


Introduction

The MyMed-Reminder simply converts a “dumb” pill box or “dumb” pill vial to a “smart pill box or “smart” pill vial. The MyMed-Reminder can remind individuals to take their medication on time, record medication history, help address the issue of forgetting or missing daily medication and notify others if the medication was taken or missed. It is suitable for home and travel use, and can be connected to smartphones or other smart electronic devices via Bluetooth. By using the MyMed Reminder app, users can access their medication records and it can also help caregivers monitor the medication situation of the user.

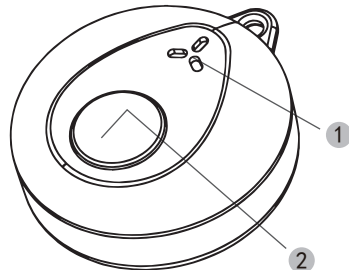

- ☆ 1 Buzzer beep sounds
- ☆ 2 Silicone button/LED light

Performance Features Include:

- ☆ Compact and portable, easy to operate
- ☆ APP with one-click location finder to prevent loss
- ☆ Automatically generate medication records
- ☆ Supports notification of guardians by email

Battery Replace

- 1 Insert a coin or another tool into the groove on the battery cover, turn counter clockwise to remove it.
- 2 Take out the old battery and insert a new button cell battery, making sure the negative terminal is facing inward.
- 3 Close the battery cover, turn clockwise with a coin or other tool, and lock the battery cover tightly. Once the battery is successfully installed, the device will automatically power on.

First Use

Before the first use, you need to open the battery cover and remove the insulation film. The factory status flashes green, indicating that the device is not paired.

Low Battery Reminder

When the battery is low, the device blinks red for 30 seconds every 10 minutes. Replace the battery.

Reminder Steps

1. Set medication information: Input the medication name, and medication time in the app.
2. Attach the device: Stick the MyMed Reminder to a commonly visible location, such as your current pill box or pill vial.
3. Connect the device: Use bluetooth to connect the MyMed Reminder to a smartphone or other smart devices.
4. Customize Reminder Alerts on the MyMed-Reminder and your phone.
5. Confirm medication: The user presses the silicone button on the MyMed-Reminder to confirm medication taken.

Medication Reminder



It's time to take the medicine!

The light flashes and the drip sounds intermittently.

Reminder Description

1. The medication schedule is transferred to the MyMed-Reminder via bluetooth and saved in its memory.
2. When it is time to take the medication, the yellow light flashes, and the buzzer beeps at the same frequency for 30 seconds. There will be 3 reminders, 10 minutes apart, for ½ an hour.
3. Press the button to stop the flashing yellow light and the buzzer reminder. The system recognizes that the drug has been taken. A green light stays on for 1 second and then turns off.
4. If the button is not pressed within the 30 minute period, the beep will stop but the yellow light will stay on for another 30 seconds. The system will record a “missed” dosage.
5. The MyMed-Reminder will transmit the “Taken or “Missed” to the APP via Bluetooth. Users can open the third-party sharing function through the APP to share the medication “Taken” or “Missed” status to a third party and to themselves by email.

Maintenance & Storage

Please keep the device in a cool, dry, and clean place, away from direct sunlight.

Specifications


Product Name: MyMed Reminder

Product Model: PR118

Front View:



Back View:

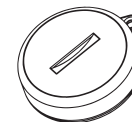

Weight: 8±2g (including battery)
Dimension: (φ30±1) × (11.5±1.0) mm
Power supply: One coin battery(CR2032)
Power Consumption: <10mA, <200uA
Working Environment:

- Temperature: 5~40℃
- Relative Humidity: 15%~93%(non-condensing)
- Atmosphere Pressure: 70~106kPa

Storage Environment:

- Temperature: -25~55℃
- Relative Humidity: ≤93%(non-condensing)
- Atmosphere Pressure: 70~106kPa

The box includes:

MyMed Reminder (including battery)X1
User ManualX1
Adhesive film.....X3

Warranty Policy

- The warranty period for this product is 12 months from the date of purchase.
- The battery is not covered under warranty.

Factory Reset

Press and hold for 3S and then quickly press 3 times to reset factory data, after reset successfully, the green, yellow, and red lights are each on for 1S and accompanied by the buzzer ringing 3 times, the factory reset will delete the original pairing information to allow the device to be paired again, and will delete all medication plans and medication records.

Non-Warranty Principles:

- Improper usage or storage;
- Damage caused by exceeding the usage conditions specified by our company;
- Malfunctions and damages caused by unauthorized disassembly or repair, and other human factors;
- Damage caused by force majeure events;
- Repairs that are not covered under warranty will incur repair and shipping costs.

FCC Declaration

Please take attention that changes or modification not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:
(1) This device may not cause harmful interference, and
(2) This device must accept any interference received, including interference that may cause undesired operation.

User Manual

MyMed Reminder

PR118

Set up APP Medication Plan

1 Download & install the MyMed Reminder app

Available in the app stores "MyMed Reminder" or scan the QR code to download and install the app, and follow the on-screen prompts to fill in the required information for registration and login in.




Android



iOS


2 MyMed Reminder connect to APP

1. Open the app and create an account or login in if you already have one account.
2. Once logged in, navigate to the device or settings section of the app.

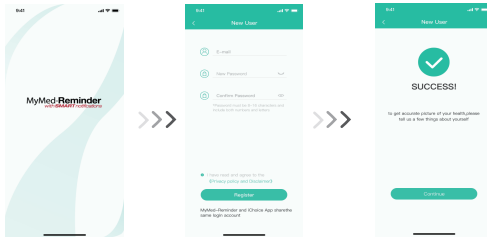


3 Login Account and Connect Device

After the login is successful, enter the device searching interface, and connect the device.
Note: Please keep your phone's Bluetooth on!

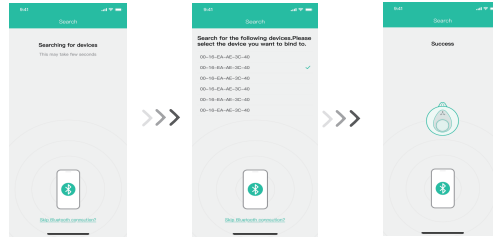


5 History Record

1. History of a medication. If you took your medication at a different time then scheduled, then you can adjust a medication record accordingly.

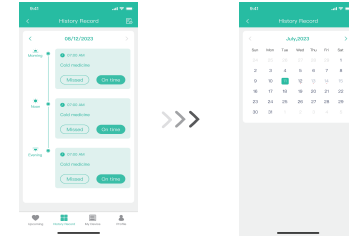


6 Set Device

- 1、Click "Customize Reminder" to display a customized notification page that can be selected in four modes.

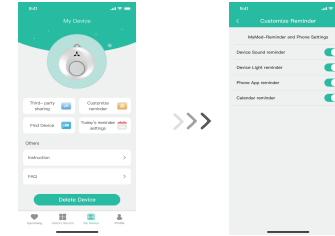


- 3、Click on "Today's Reminder Settings", choose a medication to adjust (delay or advance) only today's upcoming reminders.

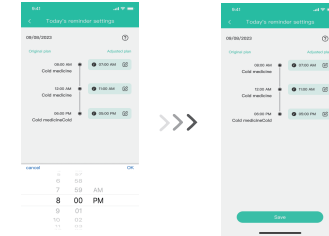


8 Profile

You can click "Profile" to change the password and other personal information.

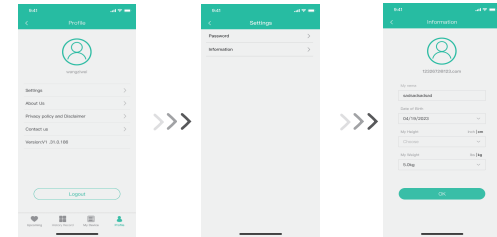

7 Third-party Sharing

Click "Third-Party Sharing", set "Email Sharing" to on, then select sharing conditions, and add an email address to share. You can enter up to xx emails. You may want to send an email to yourself.

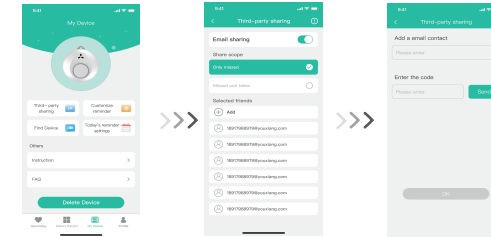

Manufacturer:
Beijing Choice Electronic Technology Co., Ltd.
2nd and 3rd Floors, No. 2 Building, No. 9 Shuangyuan Road,
Shijingshan District,
Beijing 100041, PEOPLE'S REPUBLIC OF CHINA

Issue date: September 22,2023
Version: Ver1.0