

cybersonic®

cybersonic®


Complete Sonic Oral-Care To A Higher Power.

- Cleans with gentle, yet powerful sonic strokes per minute
- Attacks heavy plaque build-up and staining
- Reaches awkward areas with ease
- Flosser attachment with pre-strung floss cartridges to make flossing easier



AMDEN CORPORATION
1111 Watson Center Rd., Suite B
Carson, CA 90745
1-800-774-7525
www.amdencorp.com

© 2005 Amden Corporation.
Cybersonic®, Cybersonic2® and Sonic White® are registered
trademarks of Amden Corporation.
40515

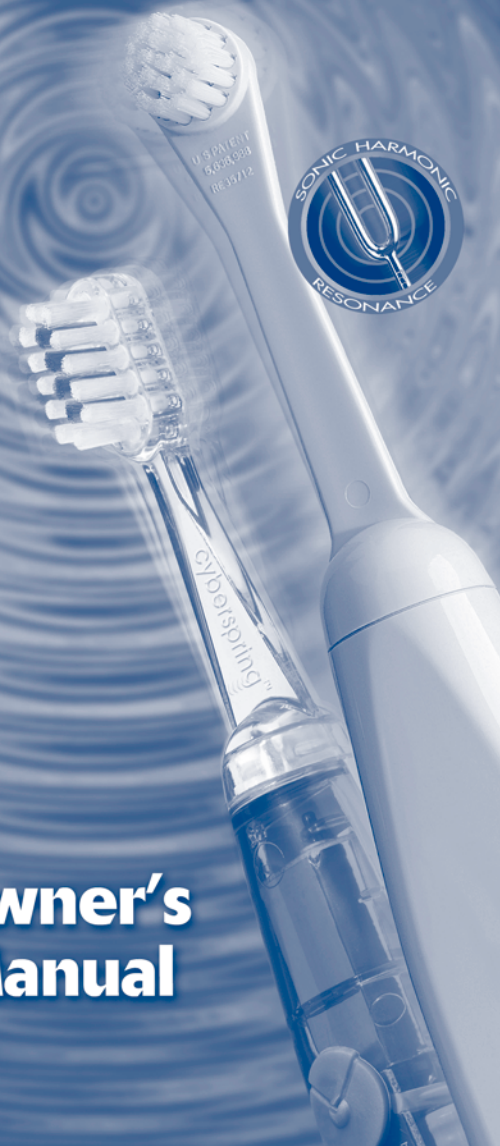

Owner's Manual

SMILE!

It Can Change Your Life.



Thank you for choosing the Cybersonic® brand, the world's first complete sonic oral-care system. Your system comes with a sonic toothbrush, sonic power flosser, sonic tongue cleaner and sonic one-step bleaching system for those confident, refreshing and dazzling-white smiles. I am proud to say that Cybersonic® and Cybersonic2® are the first and *only* "harmonic resonance" sonic oral-care systems in the world today.


It may sound complicated, but actually Cybersonic® and Cybersonic2's® technology is so advanced...it's simple. Each cleaning head is tuned, just like a musical tuning fork, to respond to sonic energy by vibrating or "resonating" back and forth, literally disintegrating plaque on contact!

Please invest the time to read these instructions and use your brush daily. If you do, I believe it will change the way you feel, the way you smile, and quite possibly, it will change your life. In order to make the most of your new Cybersonic® or Cybersonic2® system, we're ready to assist you or answer any questions toll free at: 1-800-774-7525. Or, please e-mail us at: custserv@sonictoothbrush.com

Sincerely,

Ronald K. Murayama, D.D.S.
Chairman & Founder

Your Cybersonic® or Cybersonic2® System Includes:

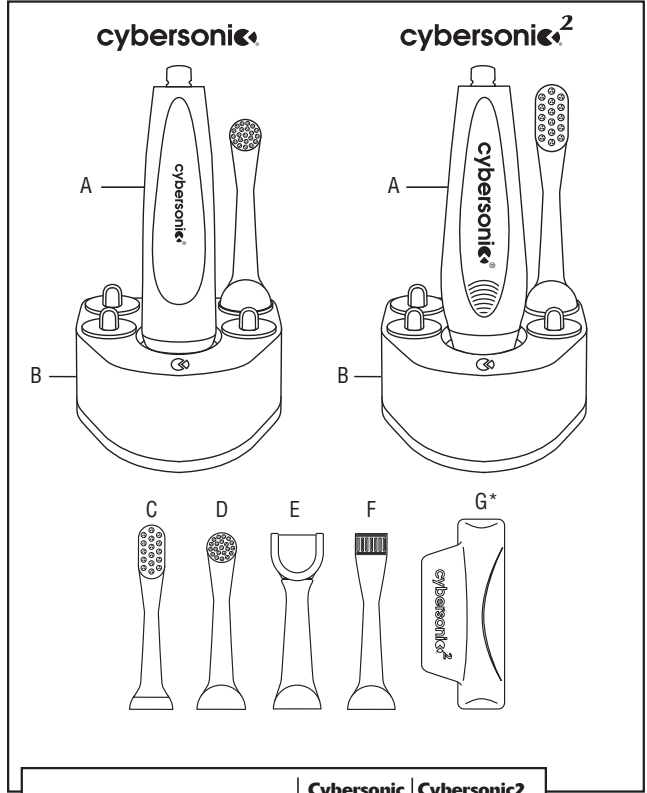

	Cybersonic	Cybersonic2
A Power Handle	✓	✓
B Charging Stand	✓	✓
C Cyberspring Brush Attachment	-*	✓
D Sonic Round Brushing Attachment	✓	-*
E Sonic Flosser Attachment	✓	✓
Pre-Strung Replacement Cartridge	✓	✓
F Sonic Tongue Cleaning Attachment	✓	✓
G Travel Case	-*	-*


Any of the attachments will fit both power handles of Cybersonic and Cybersonic2 oral care systems.

*Optional accessory.



*Case only included with C226



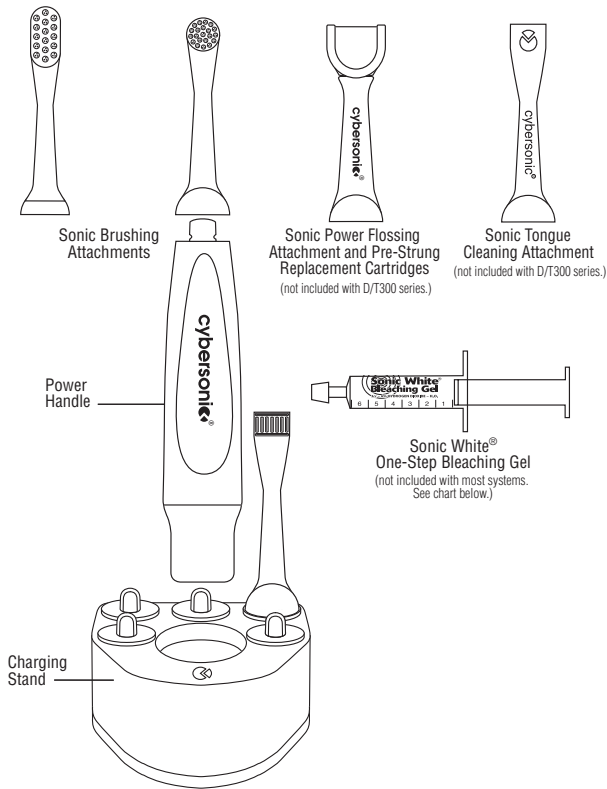

Faster Speed

- Gentle sonic strokes reduces brushing time to just one minute
- Patented sonic harmonic resonance technology
- Attacks heavy plaque build-up and staining
- Takes sonic technology to the next power



Comfortable Bristles

- Cleans hard to reach areas with ease
- Safe end rounded Dupont® Nylon bristles
- Bristle colors help identify multiple users
- Brush head cap keeps bristles clean and sanitary

Power Handle

- Comfortable to hold
- Handle won't slip out of wet hands
- Soft-touch on and off button
- Power light on Cybersonic2® indicates unit is charging

“Done-In-One™” Timer (Cybersonic2 Systems Only)


- Indicates one minute intervals with a gentle pause
- Timer starts and operates automatically when using system

Cybersonic & Cybersonic2	C201	C206	C216	C226	C266	D101	D104	D115 SERIES	DT/300
Cybersonic Power Handle	1	1	1	1	1	1	1	1	1
Charging Stand	1	1	1	1	1	1	1	1	1*
Sonic Brushing Attachment	1	1	2	2	2	1	2	2	2
Sonic Power Flossing Attachment	1	1	2	2	2	1	1	2	-
Pre-Strung Replacement Floss Cartridge	30	30	45	45	45	30	30	30	-
Sonic Tongue Cleaning Attachment	1	1	2	2	2	1	1	2	-
Sonic White One-Step Bleaching Gels	-	-	4	4	-	-	-	4	-
Q Dental Cleansing Gel	-	-	-	-	-	-	-	-	1
Decorator Cup	-	-	-	-	-	-	-	-	1
Travel Case	-	-	-	1	-	-	-	-	-

* Dome charging stand


Power Flosser

- Makes it's easy to floss every day
- Pre-strung replacement floss cartridges — simply use and discard
- Snaps on and off with ease
- Holds floss in the proper position

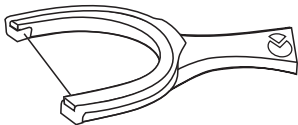


Tongue Cleaner

- Removes bad breath at the source (especially the back of the tongue)
- Acts like a mini carpet cleaner
- Flat head minimizes gag reflex
- Freshens breath naturally without mints or gum

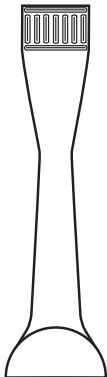


Charging Stand

- Small and compact base
- Easy to clean
- Hidden cord storage below base
- Drain hole keeps it sanitary

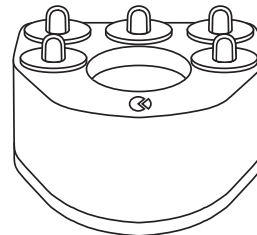

Storage Case (Optional Accessory)

- Ideal for traveling
- Stores accessories and power handle
- Quick-snap closure
- Multi-flosser compartments
- Easy to clean
- Easy to carry

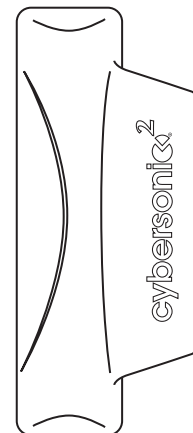

Sonic White® One-Step Bleaching Gels (Featured in select systems. See chart on page 2)

- Revolutionary one-step bleaching system
- Brushes on just like toothpaste.
- Kinetic energy stimulates the release of oxygen — the active whitening ingredient
- The easiest way to a brighter, dazzling-white smile
- Ideal for maintaining professional whitening treatments at home
- 7.5% Hydrogen Dioxide



Note: One syringe contains a one-week treatment; six days on, one day off. A one-week treatment is enough for normal yellowing or stains. Sonic White is safe for daily use.

CyberSonic® and CyberSonic2® is as easy to use as it is to maintain. Please follow these simple instructions to get started and to ensure that your system works properly and effectively for years to come. If this is your first time using a power toothbrush, you'll notice a little tingling “tickle” sensation. This is normal and minimizes after continued use.


1. Charge The Power Handle 24 Hours.

Before using your system for the first time, it is very important to charge the power handle (for at least 24 hours before for the first use). Failure to do so could damage the battery system.

- Plug charger into outlet
 - Place power handle into the stand, the charging indicator will illuminate on CyberSonic2® units.
 - Charging stand and power handle may warm slightly
- After the initial charging period, you may begin using your system.

2. Install or Remove Attachment.


Select a cleaning attachment. Push straight onto the end of the power handle until you hear a click. For optimum results, the face of the brush or cleaning attachment should be on the same side as the ON/OFF button. However, you may orient the attachment position anyway you feel comfortable. Do not try to pull the attachment straight off. To remove, push the attachment to the side while pulling at the same time.

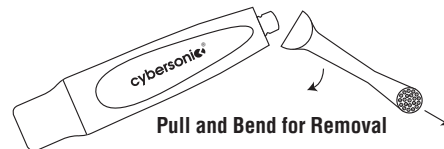

3. Power On.

Locate the push button on the side of the power handle. Push this button to power on and again to power off when you're finished cleaning.

The power turns off automatically when the power handle is placed into the charging stand. The handle is designed so you can use it anyway that feels comfortable in your hand.

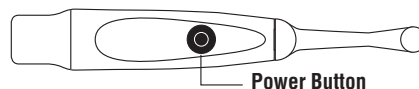

Charging the Batteries

Keep the charging stand plugged in continuously. And for optimum performance you should place the power handle back into the charging stand after every use. You cannot overcharge the batteries. To extend battery life, unplug the charging stand every other month. Plug the charger stand back in after 48 hours or if the handle runs out of power.

Batteries should last three to five years depending on how frequently you use the product.

Battery Disposal

We use a nickel-cadmium rechargeable battery. Please do not dispose of this appliance in the household waste basket. Instead, dispose at a recycling center.

Cleaning the System

After every use, rinse and clean any toothpaste or residue from the attachments and the easy-grip power handle. Take off the attachments and rinse parts separately, then wipe them dry.

Replacing the Brush Attachments


We recommend replacing the brushing attachments every three months. Replacement brushes may be purchased at your local retailer or by calling: 1-800-774-7525.

Using the System



The System is an advanced sonic oral-care system that disintegrates plaque on contact and features a full range of attachments for complete oral health. Follow the simple instructions below to begin using the attachments properly.

Sonic Brushing Attachments

1. Place the Sonic brushing attachment onto the power handle. Remove the brush head cap and apply toothpaste or gel if desired.

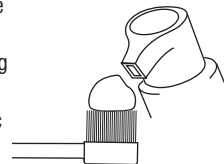


2. Using moderately firm pressure, place the bristles on your teeth. Turn on the easy-grip power handle. Move the brush head evenly over all surfaces of your teeth. No extra scrubbing or brushing motion is required in brushing your teeth and gums. Just press firmly and the sonic technology does the rest.
3. When brushing the sides of your teeth, place the brush head at a 45° angle. Apply firm pressure while moving along the gum line. The "Hyper-Bass" technology directs sonic cleaning below the gum line. Spend at least one minute brushing all surfaces of your teeth.

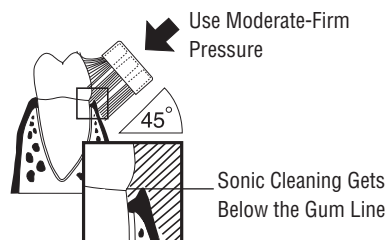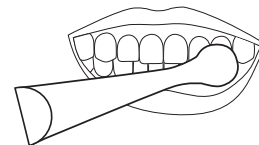

Tongue Cleaning Attachments

1. Place the sonic tongue cleaning attachment onto the power handle and use with or without paste/gel.

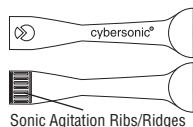

2. Place the tongue cleaner into your mouth and on top of your tongue. Use light pressure and rub over the top surface in back and forth motions. Pull or scrape the residue out and expel.



Sonic White® One-Step Bleaching Gels

1. Apply the gel onto the brush head. Push the plunger in until it reaches the first dosage marker. Each notch equals one day's treatment (replace end cap on syringe).

2. Spread an even film of gel across your teeth and work the gel using small circular motions.



Continue for one to two minutes for stubborn stains. Rinse and look for results within one to three weeks.

*Contains no animal by-products. Deionized water, Hydrogen Dioxide 7.5%, Carbopol and Triethanolamine.
 Note: One syringe contains a one-week treatment; six days on, one day off. A one-week treatment is enough for normal yellowing or stains. Continue for additional one to two weeks for more stubborn stains. Sonic White is safe for daily use. For best results, do not bleach teeth right after using a fluoride toothpaste or rinse. During normal use, minor gum or tooth sensitivity may occur or white film or spots may appear, but are harmless and temporary. If you experience prolonged sensitivity, pain or bleeding gums, stop use and consult your dentist.

Sonic Flosser Attachments Plus Pre-Strung Replacement Heads

1. Place the flosser attachment onto the power handle. Then snap the pre-strung replacement head into the flosser attachment.

2. Turn on the power handle. Insert floss between teeth. For tight teeth, you may need to move the flosser from side to side.



3. To clean between your teeth, slide the floss up and down both sides of your teeth. Also floss under the gum line.



4. To replace the used flossing head, place thumbnail under base and snap upward and out.

5. Dispose of replacement heads after each use. Order replacements by calling 1-800-774-7525.

Commonly Asked Questions

Q. Will the system loosen crowns, fillings or other permanent dental work?

A. No, they are designed to be used for all dental work.

Q. How often should I replace my brush heads?

A. Dentists recommend changing every three months. Your brush may not look worn but it is, and a worn brush does not clean well.

Q. Is it safe to use for a person with a pacemaker?

A. It is always advisable that a person with a special medical conditions, including pacemakers, consult with a physician or pacemaker manufacturer.

Q. Can I travel with my power toothbrush?

A. Yes. A fully charged brush will hold a charge for at least a week of normal use, and the brush head cap keeps the bristles clean.

Q. Is it suitable for kids?

A. Yes. The original round brush is actually a perfect size for children. You may start your kids using the system at the same age you would a manual toothbrush. Parental supervision is always advised with a new oral cleaning product.

Q. Why do I sometimes need to push the ON/OFF button more than once?

A. The electronic switch takes 1 second to re-set before it will respond to another push. Simply wait 1 second before pushing the button again.

Q. How do I purchase replacement attachments?

A. Replacement attachments may be purchased at your local retailer or by calling 1-800-774-7525. You may also access our web site at www.amdencorp.com.

Important Safeguard Precautions!

Please Read Before Use.

Our toothbrushes have been carefully engineered for effectiveness and safety. However, basic safety precautions should always be followed when using any electrical product.

To reduce the risk of electrical shock, fire or injury:

- Avoid placing the charging stand where it might fall or be pulled into a tub, shower, sink, or toilet.
- Do not use the system while bathing or showering.
- Do not reach for charger if submerged in water. Unplug the stand immediately and drain the water source before retrieving it by the cord.
- Keep the cord away from heated surfaces.
- Do not use the system if the charging stand, cord, or plug are damaged or have been immersed in water. Call Customer Service at 1-800-774-7525 for further instructions.
- Supervise children under 12 and disabled people when they are using this electrical device.

Warnings

- Do not use toothpaste that contains baking soda as they will react with the ABS plastics that your system is made from, causing cracks or breakage.
- Consult your dentist before using product if you've had periodontal disease, dental abscesses, oral infections, oral pain, oral discomfort, bleeding gums, loose teeth or toothaches.
- Consult your physician before use if you are a cardiac patient with a pacemaker and have concerns. However, the system produces only a small magnetic field in the power handle — no stronger than emitted by an electric shaver.
- If you experience persistent pain, bleeding or discomfort during or after use, consult your dentist.
- Not intended to replace regular dental check-ups or cleanings.

One Year Limited Warranty

(excludes toothbrush bristles and replacement floss cartridges)

The Cybersonic® and Cybersonic2® systems Limited Warranty covers mechanical and/or electrical failures of the power system. Under warranty terms and conditions, the system shall be repaired or replaced at the manufacturer's discretion. Defective system and proof of purchase must be submitted along with warranty claims. This warranty does not cover:

- Brush heads
- Floss cartridges
- Damage caused by misuse, abuse, neglect or alterations
- Normal wear and tear, including chips, scratches, abrasions, discolorations or fading

Warranty claims and proof of purchase should be forwarded to:

Amden Corporation

1111 Watson Center Rd., Suite B

Carson, CA 90745

1-800-774-7525

www.amdencorp.com

