



Step 4. Connecting your Home Phone

Connect your home phone cord to the gray connection on the back of the Console.

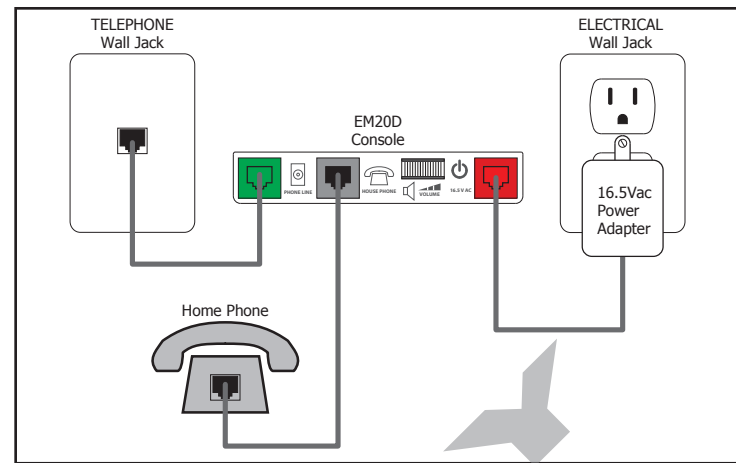


Step 5. Setting up your Help Button

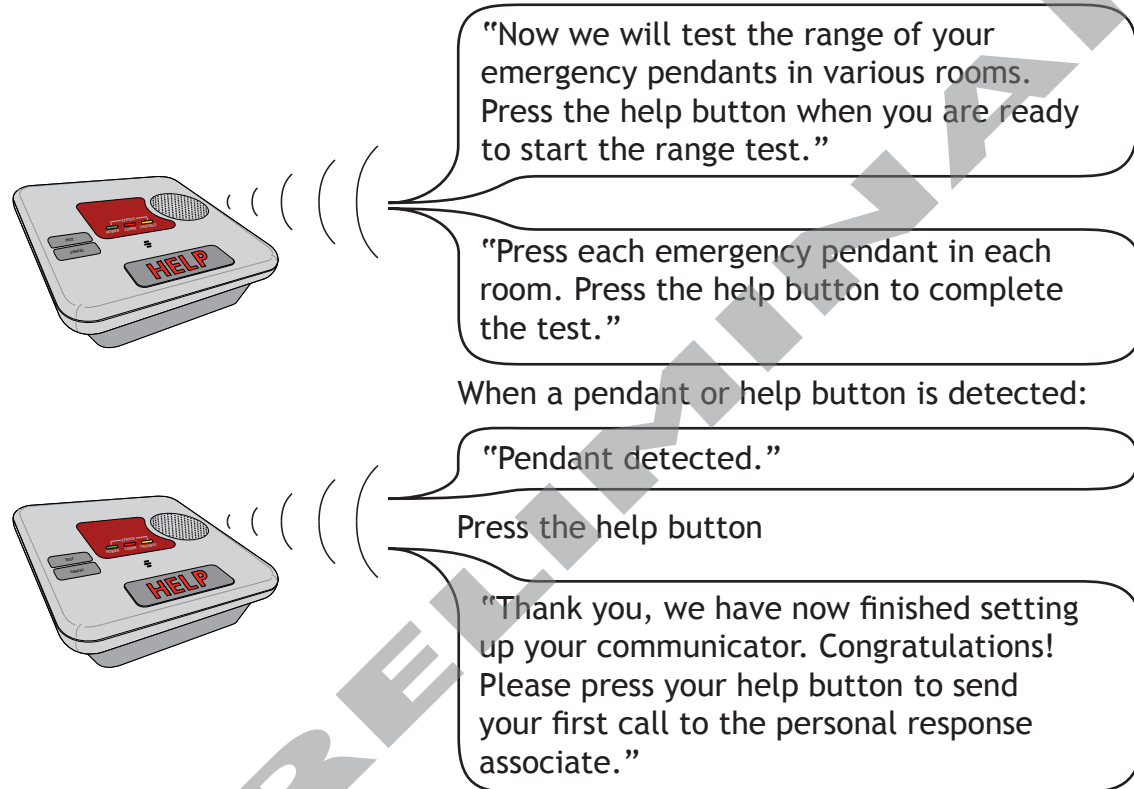

Step 6. Placing a Test Call

Press and release the Test Button within 3 seconds to place a test call.

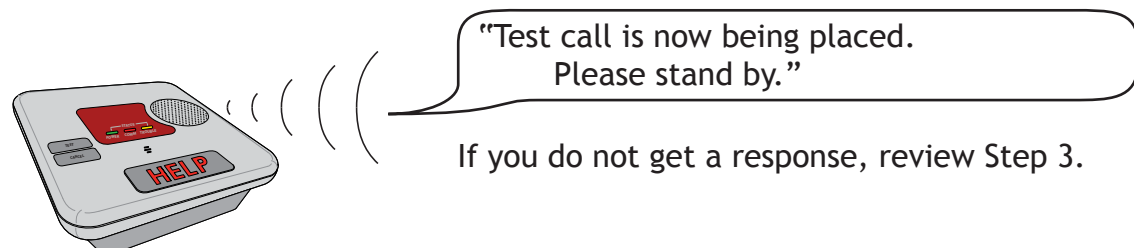


LifeStation[®] Quick Start Guide

SELECTING YOUR HELP BUTTON OPTION | CONNECTING THE POWER SUPPLY | CONNECTING THE PHONE LINE | CONNECTING YOUR HOME PHONE | SETTING UP YOUR HELP BUTTON | PLACING A TEST CALL



What is Included

- One EM20 Console
- One 12V DC Battery - installed
- Color-coded phone cord
- 16.5VAC Power Adapter
- 1148 Emergency Pendant Wrist strap
- Break-away Lanyard

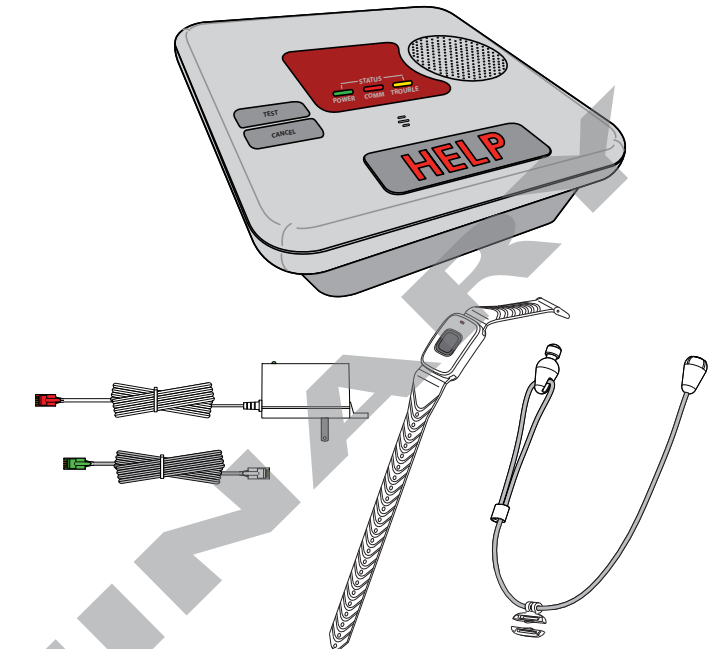

Step 1. Selecting your Help Button Option

Choose either the wrist strap or break-away lanyard option for the 1148 Emergency Pendant.

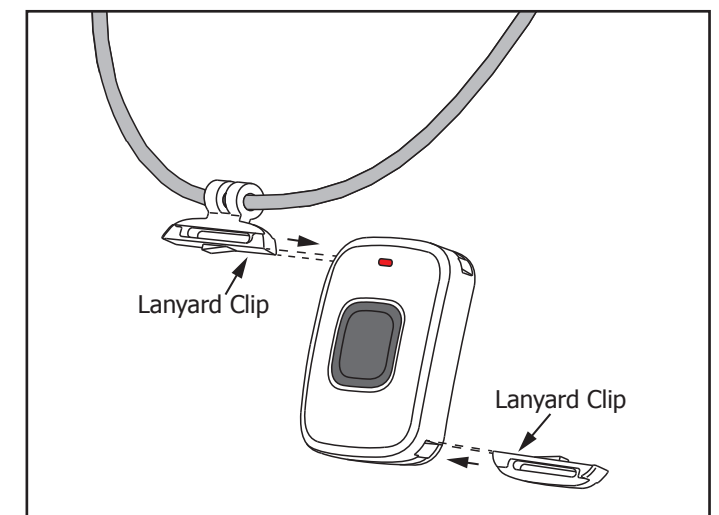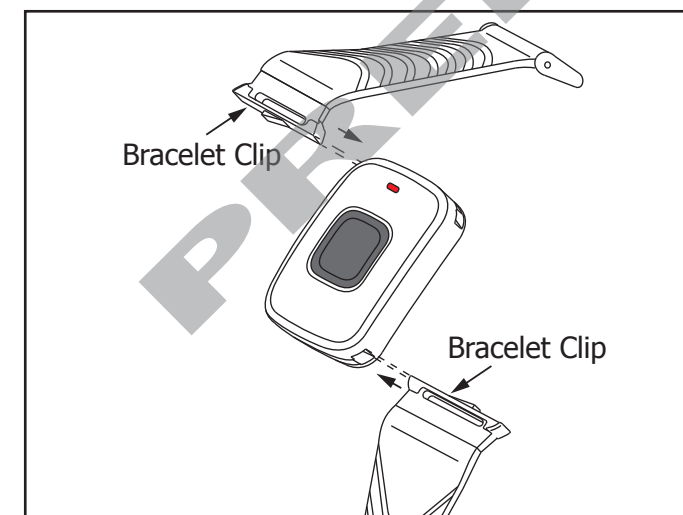

LT-1329 © 2014 LifeStation 14475

Step 2. Connecting the Power Supply




The Voice Prompt Setup begins automatically when power is applied. The setup cannot be exited until completed. If power is somehow interrupted and the setup does not complete, it will start over when power is restored. During the Voice-Prompt Setup the last message will repeat once for each minute of inactivity.

1. Unpack the power supply and fully unroll the cord.

2. Connect the power supply to the red connection on the back of the Console.

3. Connect the power supply to an electrical outlet and secure with the provided screw through tab on the power supply.

When power is applied, the Voice Prompt Setup begins.

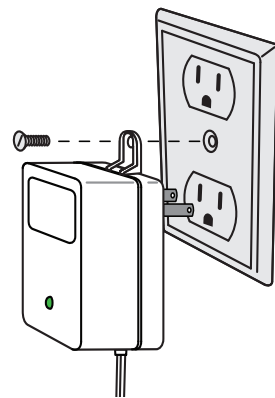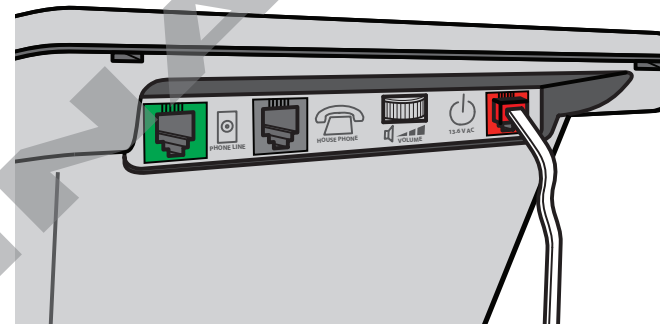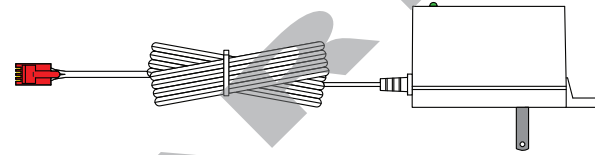

"Welcome! We will now help you setup your communicator step-by-step. If you miss a step, wait, and we will repeat the instructions."




"Please find the quick start guide located in the box. When you have the guide in front of you, press the help button."

Step 3. Connecting the Phone Line

1. Unpack the color-coded phone cord and fully unroll.

2. Connect the color-coded phone cord to the green connection on the back of the Console.

3. Connect the opposite end to the phone wall jack.

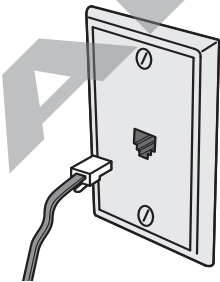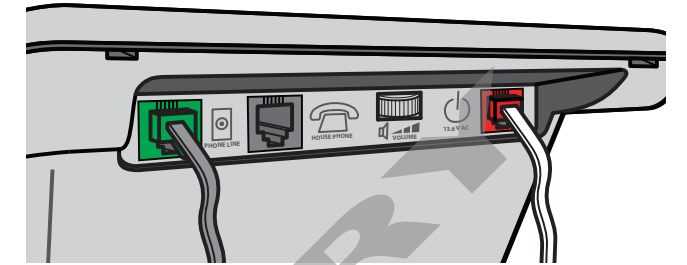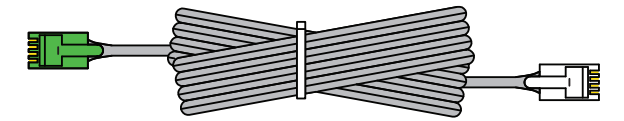

Next, we will plug the gray phone cord into a phone jack in the wall. You may need to unplug your phone temporarily, and reconnect it in the next step."

If no phone line is detected within one minute:

"Your communicator has not been correctly plugged in to the phone line. See step 3 in the quick start guide."

When a phone line is detected:

"Good, your communicator is now plugged in. You may press the help button to continue."



Next, we will plug your phone into the back of the communicator if necessary. Refer to step 4 in the quick start guide. Please press the help button to continue."