Live!y™

Live!y™

www.mylively.com

© 2014 Lively. All rights reserved. Lively and LivelyGram are a trademark of Lively, Inc.

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CONTENTS (english language)

- p. 2 Safety watch overview
- p. 4 Installing the battery
- p. 5 Activating the watch
- p. 6 Pushing the help button
- p. 8 Receiving medication reminders
- p. 10 Having fun with step-counting
- p. 11 Using the safety watch with a smartphone
- p. 12 Customizing your safety watch
- p. 13 Using Lively's custom watch band

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IMPORTANT

You must activate your Lively system, at app.mylively.com/start, in order to set up the safety watch. During the activation process, you will receive a test emergency phone call. Once activated, your safety watch must maintain a constant connection to the Lively hub at home, or be paired on-the-go with a smartphone, for emergency-response use.

Welcome to emergency response, *reimagined.*

Lively is much more than a traditional emergency response “button.” Of course, Lively is there when you need to summon help at the push of the button. But we wanted your safety watch to be useful—and fun—for all aspects of your active lifestyle.

We started with a fashionable watch that you can customize and accessorize to fit your personality. With smart medication reminders, the safety watch notifies you only when you might have missed taking medication or vitamins. And the watch includes a pedometer for simple step-counting to help you measure how active you are being throughout the day.

Enjoy and stay Lively.

- **Waterproof**
You can use the Lively safety watch in the bath and shower.

- **Help button**
Press to receive a call from the Lively Care Team to the phone number you designate.

- **Day & night display**
A large, crisp display with back lighting for easy viewing.


 **Tip:** Tap the edge of the watch face to light up the display.

- **Interchangeable bands**
Your Lively safety watch comes with an easy-to-use, comfortable band, which can also be replaced with any standard watch band (20mm width).

Installing the battery

- Remove the Lively safety watch and battery from its packaging.
- Insert the coin battery (standard Li-Ion battery, CR-2032) into the watch.
- Using a coin, twist the back cover firmly to close it. *This is important to ensure waterproofing.*

 **Tip:** After a new battery has been installed, the safety watch will re-activate once it is within range of the Lively hub.

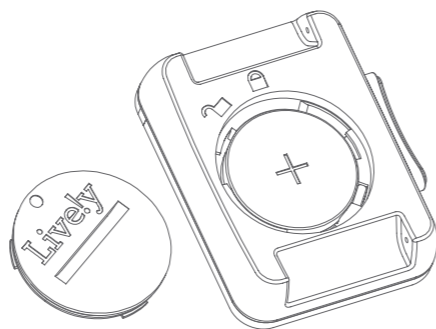Activating the watch

- Make sure the safety watch is near the Lively hub for activation.
- Make sure you have completed the online Lively activation process up to the “Activate system” step.

- Follow the on-screen instructions to activate the watch.

During this process, you will perform a test of the help button and receive a phone call from the Lively Care Team.

If you do not currently have online access via a computer or tablet, please call the Lively Care Team at 888-757-0711.

- After you have completed activation and successfully received a test phone call, your watch will be ready for use!

Note: The Lively safety watch will not function until you complete the activation process.

 **Tip:** At a minimum, we recommend testing your help button when you replace the battery. Just press the button to test the response of Lively’s Care Team. Once you receive a confirmation phone call in response, let Lively know you’re ok and there’s no need for assistance.

Pushing the help button

Important: Lively is not a substitute for 9-1-1. In the event of a critical emergency, always contact 9-1-1 immediately. By subscribing to the Lively service, you accept, without limitation or qualification, the Lively Terms and Conditions.

- Your safety watch must be connected to the Lively hub or paired on-the-go with your smartphone for emergency response use. See page 12 to pair with a smartphone.
- If you need assistance, press the help button.
- The Lively Care Team will call you at the phone number you provided during the initial setup of Lively. You can update this phone number by visiting app.mylively.com from a computer or tablet.


Tip: The watch's range of use in the home or the yard is up to 1,500 ft., when unobstructed.

- The watch will display the following while the phone call is being placed to you. >>
- If you do not answer the telephone, or if you request emergency services while talking to a Lively Care Team member, emergency services will be dispatched.

Note: Emergency services will be dispatched to the address you provided during the initial setup of Lively. To update your home address, visit app.mylively.com

- The watch will display the following while help is on the way. >>

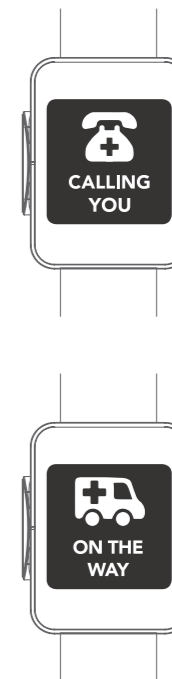Receiving medication reminders

- With Lively's safety watch and pillbox activity sensors, you will be reminded automatically whenever you've missed taking your medication.
- Set up medication notifications during the activation process, or go to app.mylively.com and choose Settings > Set Notifications. Simply choose when you want to take medication, and start receiving notifications.
- If the pillbox sensors do not detect movement during the window of time when the medicine should be taken, our Lively system automatically sends a helpful reminder to the watch. >>

- The medication reminder will repeat once after fifteen minutes, unless the pillbox sensor has registered movement, indicating the medication has been taken.
- If medication does not register as taken by the end of the allotted time, notifications will be sent to the Circle members who have been set up to receive them.

Having fun with step-counting

- Your Lively safety watch also includes a pedometer to keep track of your steps, wherever you are, throughout the day.
- Activity is sent to the Lively hub approximately every 15 minutes, when you are in the home, and when you return from being away.
- To see how many steps you take each day, simply tap on the edge of the watch face to see your latest activity.
- You and your family can see step activity by using a Lively smartphone app or visiting app.mylively.com on a computer or tablet.

Using the safety watch with a smartphone

- Connecting your Lively safety watch to a smartphone enables 24/7 help outside the home.
- To use the watch with a smartphone, download the Lively app and sign in. Smartphone use requires Bluetooth LE and GPS; please go to support.mylively.com and search “smartphones” to find a list of supported phones.
In order for the alert button to work, your watch must be within range of the phone. Your watch screen will confirm connection. >>
- 911 services will be dispatched to the GPS location automatically provided by your phone, unless you confirm help is not needed.

For more information on how to operate your safety watch with a smartphone, please visit support.mylively.com


Customizing your Lively safety watch

The Lively safety watch can be customized by logging in at app.mylively.com from a computer or tablet. Choose Settings > Safety Watch.

In settings you can:

- Pick analog or digital watch display
- Turn on/off vibration of the watch
- Turn on/off medication reminders

Visit Settings at app.mylively.com to see more customization options.

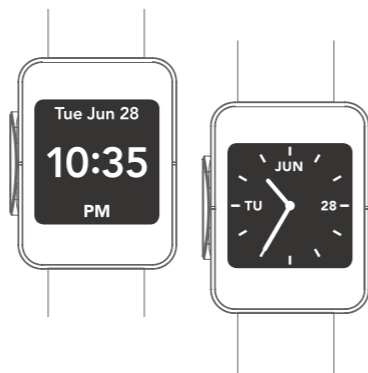Using Lively's custom watch band

We hope you enjoy the elegant custom watch band from Lively. We've taken great care to design a watch that provides the ultimate comfort, convenience and durability.

Features of our custom band include:

- Wide-body bands made of super soft-touch silicone
- Waterproof and washable
- Oversized holes for easy adjustments

Feel free to replace the Lively band with your own favorite band, as long as it has a standard width (20mm).

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