READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT

Keep this manual in a safe location for future reference

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CONTENT

1. IMPORTANT SAFETY INSTRUCTIONS.....	2	8. Evaluating Results	5
2. Introduction	3	9. Other functions.....	5
3. Package Contents	4	10. Care and Maintenance	6
4. Parts and Controls.....	4	11. Troubleshooting Guide	6
5. Battery Installation	4	12. FCC Compliance Information	7
6. Initial Set Up for Automatic User Recognition with App.....	4	13. Technical Specifications	8
7. Taking Measurements	5	14. Warranty	8

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact customer service. Please see warranty for service contact.

Dear Customer,

Thank you for choosing one of our products. Please read these instructions for use carefully and keep them for later use. Make them accessible to other users and observe the information they contain.

1. IMPORTANT SAFETY INSTRUCTIONS

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ AND UNDERSTAND THIS ENTIRE MANUAL INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.

WARNING

WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury.

CAUTION

CAUTION indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.

NOTICE

NOTICE addresses practices not related to personal injury, such as product and/or property damage.

WARNING – To reduce the risk of fire, electric shock, or serious personal injury:

- To prevent falls, do not step onto the scale with wet feet or if the surface of the scale is damp.
- Do not stand on the outer edge of the scale, otherwise it may tip.
- Keep children away from packaging materials (risk of suffocation).

CAUTION – To reduce the risk of personal injury or product/property damage:

- This product is intended solely for household use and not for medical or commercial purposes.
- Protect the device from impacts, moisture, dust, chemicals, large temperature fluctuations, and nearby sources of heat (ovens, heaters).
- Do not use on children under 10 years of age.

NOTICE

- Place the scale on a hard, flat, level surface; a hard surface is required for accurate measurements.
- Measurement variations are normal, as this scale is not calibrated for use in a professional medical environment.
- The device should be cleaned from time to time. Do not use any abrasive cleaning products and never submerge the device in water.
- Make sure that no liquids come into contact with the scale. Never submerge the scale in water. Never rinse it in running water.
- Do not place any objects on the scale when not in use.
- Do not press the buttons violently or with pointed objects.
- Do not expose the scale to high temperatures or strong electromagnetic fields.
- The scale's capacity is 400 lb (181.4 kg).
- The scale is initially set to weigh and measure in 'lb'.
- Repairs must only be performed by Customer Service.
- The accuracy of the measurements and service life of the device depend on its careful handling.

Battery Handling Safety Precautions

- Use only the size and type of batteries specified.
- Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- If the battery in the device is depleted or the device will not be used for a long period of time, remove the battery to prevent damage or injury from possible battery leakage.
- Do not try to recharge batteries not intended to be recharged; they can overheat and rupture (follow battery manufacturer's directions).
- Do not dispose of batteries in fire, batteries may explode or leak.
- Clean the battery contacts and also those of the device prior to battery installation.
- Remove discharged batteries from the product and dispose/recycle in compliance with all applicable laws.
- Keep batteries away from children and pets. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately.

2. Introduction

This *Bluetooth*[®] glass scale is intended for private use at home only.

The scale has the following functions that can be used by up to eight persons:

- Body weight
- Body Mass Index (BMI)

The scale also includes the following features:

- Switches between pounds (lb) and kilograms (kg)
- Automatic shutoff
- Low battery indicator
- Automatic user recognition
- Stores last 30 measurements for each of eight users
- The scale uses *Bluetooth*[®] low energy technology and transmits via the 2400 - 2480 MHz frequency band, transmission power max. 0 dBm.

System requirements for the 'HealthCoach' app

Bluetooth® ≥ 4.0, iOS from version 8.0, Android™ from version 4.4 with *Bluetooth*® low energy technology.

3. Package Contents

- 1 x GS 435 *Bluetooth*® Glass Scale
- 1 x 3V CR 2032 battery
- 1 x Instruction manual

4. Parts and Controls

Front:


- 1. Display

Top:

- 2. Up/Down Buttons


Rear:

- 3. Unit /Reset Button
- 4. Battery Compartment

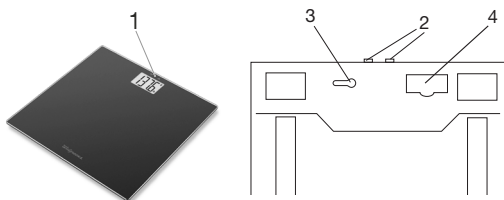Display:

- 5. *Bluetooth*® symbol for connection between scale and smartphone
- 6. Scale User
- 7. BMI
- 8. Body Weight Unit
- 9. Body Weight

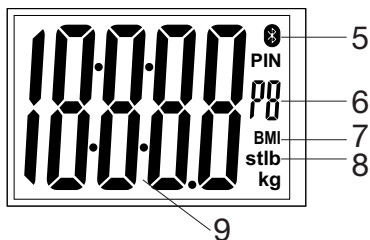5. Battery Installation

Remove the Battery Compartment Cover and insert one 3V CR 2032 battery (included) according to the indicated polarity. If the scale does not function, remove the battery and reinsert it, ensuring proper polarity. Replace Battery Compartment Cover. "APP" will flash on the display if the scale memory is empty.

6. Initial Set Up for Automatic User Recognition with App

For a correct automatic user recognition, you must save the personal user data on the scale via the app. The scale has 8 user memories which allow, for example you and your family members to save personal settings. Place the scale on a hard, even floor; a firm floor covering is required for correct measurement.

- Activate *Bluetooth*® in the smartphone settings. Keep the smartphone close to the scale in order to maintain an active *Bluetooth*® connection during initial set up.
- Install "HealthCoach" from Apple App Store or Google Play.
- Start the app and follow the instructions in the app to add your scale to the app.

NOTE: To select the GS 435 in the app, the scale must be on for each synchronization. When the scale is off, switch on the scale by briefly stepping onto the platform of the scale.

- For automatic personal recognition, the first measurement must be assigned to your personal user data. Follow the app instructions for this. Step onto the scale and ensure that you are standing still with equal weight distribution and with both legs on the standing surface.

7. Taking Measurements

- Place the scale on a hard, flat surface with no carpet or rugs; a hard surface is required for accurate measurements.
- Step on the scale. Stand still with your weight distributed equally on both feet. The scale begins to take measurements straight away.
- The weight is displayed.
- If a user has been assigned, the correspondent user (P1 to P8) and the BMI are displayed. This measurement is saved for *Bluetooth*[®] synchronization.
- If the user is not recognized, only the weight is shown. Press the “Up / Down” buttons on the top of the scale to select a user manually. Step on the scale again.

Tips

- If possible, always weigh yourself at the same time every day (ideally in the morning), after having used the restroom, on an empty stomach and with as little clothing as possible, in order to obtain comparable, consistent results.
- Stand up straight and still while weighing.

Switching off the scale

After displaying your weight, the scale will switch itself off automatically.

8. Evaluating Results

Body mass index (BMI)

The body mass index (BMI) is a number that is often called upon to evaluate body weight. The number is calculated from body weight and height, the formula is – body mass index = body weight : height². The measurement unit for BMI is [kg/m²]. According to the BMI, weight is classified for adults (20 years and over) using the following values (source: WHO):

Category		BMI
Underweight	Severely underweight	< 16
	Moderately underweight	16–16.9
	Slightly underweight	17–18.4
Normal weight		18.5–24.9
Overweight	Pre-obese	25.0–29.9
Obese (overweight)	Class I obese	30–34.9
	Class II obese	35–39.9
	Class III obese	≥ 40

9. Other functions

User assignment

It is possible to assign up to eight users. For a new measurement, the scale assigns the measurement to the user whose most recently saved measurement falls within +/- 6.6 lb (3 kg) of this.

Saving the assigned measurements on the scale

A total of 30 measurements per user can be saved on the scale when the app is not active. The saved user measurements can be transferred to the app if you switch on the scale and open the app within *Bluetooth*[®] range. It is not possible to transfer data when the scale is switched off.

Add Existing Scale User after Reinstallation of app

If a scale user already exists (e.g. from an old app version) this scale user can be integrated into the new reinstalled app by inserting the 4-digit user PIN into the app. The app is asking for this user PIN when required. The PIN is displayed on the scale's LCD display when the app is connected. Press the "Up / Down" buttons to select a user and wait 3 seconds. Briefly press the Unit/Reset button.

If the integration was not successful, delete the scale user and create a new scale user.

Deleting scale data

To delete scale user data including measurements, switch on the scale, select the user manually by pressing the "Up / Down" button and press and hold the "Unit / Reset" button for five seconds. 'dEL' will appear in the LCD Display.

To delete all measurements and user data on the scale, switch on the scale and press and hold the "Unit / Reset" button for five seconds; 'CLr' will appear in the LCD Display. You then have to carry out the steps in chapter 6 "Initial Set Up for Automatic User Recognition with App" again.

Use scale without user measurement

If the user does not want to initialize the scale using the app or does not want to use the app at all, the user can stop "APP" flashing by shortly pressing the "Unit/Reset" button. Afterwards, the scale can be used as a normal weight scale without any user memory function.

10. Care and Maintenance

Replacing batteries

Your scale is equipped with a low battery indicator. If the battery is weak, 'Lo' will appear in the LCD Display and the scale will automatically shut off. In this case, replace the battery with a 3V CR 2032 battery.

Cleaning

Clean the scale periodically using a damp cloth and, if necessary, a small amount of mild detergent.

Disposal

Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

11. Troubleshooting Guide

If the scale encounters an error during measurement, the following is displayed:

Display	Cause	Remedy
Only weight is displayed.	Unknown measurement as it is outside the user assignment limit or a unique assignment is not possible.	Press the 'Up' button and wait three seconds to choose a user manually. Repeat measurement. Use the app if no user exists.
Lo	The maximum load-bearing capacity of 400 lb (181.4 kg) was exceeded.	The load must not exceed 400 lb (181.4 kg).
Incorrect weight is displayed.	No flat, stable surface.	Place the scale on a flat, stable surface.
Incorrect weight is displayed.	Not standing still.	Stand as still as possible.
Incorrect weight is displayed.	Scale zero setting is incorrect.	Wait until the scale switches itself off. Activate scale, wait for '0.0 lb' to appear and repeat measurement.
"APP" is flashing.	The scale needs to be connected with the app.	Connect with app.

No <i>Bluetooth</i> ® connection (Ⓜ symbol is missing).	Device outside the range.	The range is approx. 32,8 ft (10 m). Walls and ceilings reduce the range. Other radio waves may interfere with the transmission. Therefore, do not position the scale near devices such as WIFI routers, microwaves, inductions hobs etc.
	No connection to the app.	<ul style="list-style-type: none"> • Briefly step on the scale platform to switch on the scale and open the app. • Close all apps completely (including in the background) and reopen the app. • Switch <i>Bluetooth</i>® off and back again. • Switch the Smartphone off and back again. • In the <i>Bluetooth</i>® settings of the operating system, the scale must not be connected. • Remove the batteries from the scale for a short period and reinsert them.
FULL	User memory space is full. No more measurements can be saved.	Open the app and switch on the scale. The data is transferred automatically. This may take up to one minute.
Lo	The battery in the scale is empty.	Replace the battery of scale.

12. FCC Compliance Information

GS 435 *Bluetooth*® Glass Scale

Responsible Party – U.S. Contact Information
 Beurer North America LP
 900 N Federal Highway, Suite 300
 Hallandale Beach, FL 33009
 United States
 1-800-536-0366

FCC Compliance Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:
 1. This device may not cause harmful interference, and 2. This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

13. Technical Specifications

Battery type	1 x 3V CR 2032 battery
Measuring range	6.6 lb to 400 lb (3 kg to 181.4 kg)
Weight unit / increments	lb or kg / 0.2 lb increments (100 g)
Product Dimensions	11.81 x 11.81 x 0.83 in (300 x 300 x 21 mm)
Product weight	3.2 lb (1450 g)
Subject to technical modifications.	

14. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Walgreens *Bluetooth*[®] scale, model GS 435, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Walgreens *Bluetooth*[®] scale, model GS 435, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at **1-800-536-0366** or at **info@beurer.com** to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

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Fabriqué en Chine

The *Bluetooth*[®] word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Beurer North America is under license. Other trademarks and trade names are those of their respective owners.

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Google Play and the Google Play logo are trademarks of Google LLC.

Android is a trademark of Google LLC.

Qualified Design ID (QDID): 100518