




**One To One™ Series  
Caregiver Alert AT-500  
Instruction Manual**

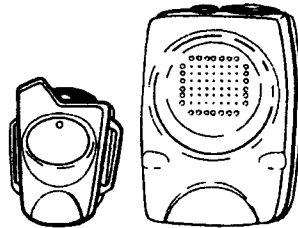

(Actual style of transmitter and receiver may vary from illustrations):


**Your Caregiver Alert includes**

- Battery operated receiver w/ vibration and audio alerts
- Waterproof transmitter, battery included
- Wrist strap & Pendant strap for transmitter
- Belt clip for receiver


**Installing batteries in receiver:**

You'll need to buy 3 "AA" batteries for the receiver.  
In typical use alkaline batteries will last up to 6 months.


1. To open the battery compartment, first remove the belt clip by sliding clip down.

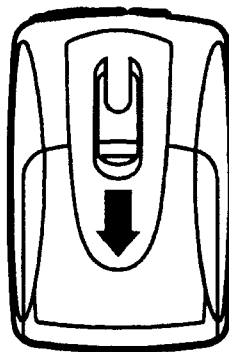


2. Open battery compartment by sliding door in direction of arrow on door.

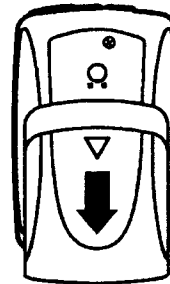


3. Install the 3 alkaline "AA" batteries making sure the batteries are aligned properly.

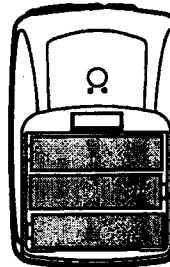

4. Slide door onto receiver making sure the door snaps into place.

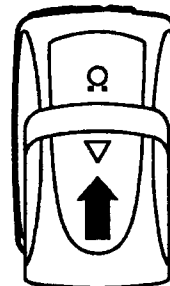


5. Slide Belt clip onto back of receiver; note that the clip can be snapped into place in 4 directions 90° of each other.

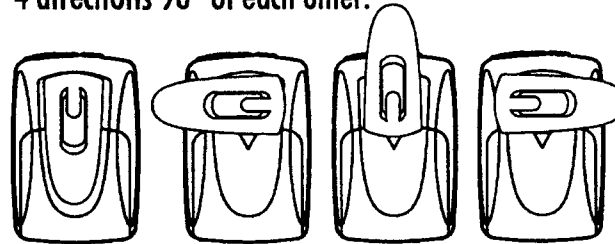


**Replacing battery in transmitter:**

(Note transmitter comes with Lithium battery- LR2032 installed)


1. To open battery compartment door, remove screws in door using a small Philips head screwdriver.

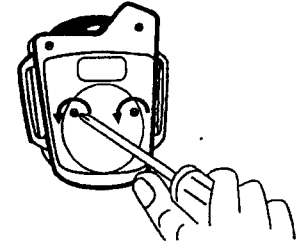


2. Carefully slide old battery (note battery alignment) out and discard as per local codes.

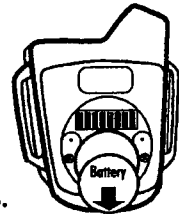

3. Slide new Lithium battery in place in same alignment as old battery.

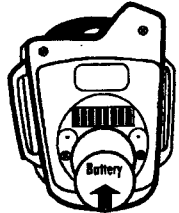


4. Place battery compartment door on to transmitter and screw in place.

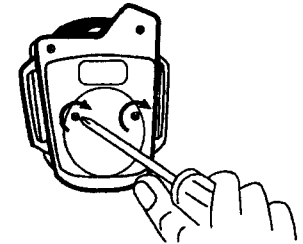

Test Caregiver Alert by pressing button on Transmitter to verify the transmitter and receiver function properly.

### Operation

Your One to One™ Caregiver Alert receiver has the options of vibration, sound or both vibration and sound by simply sliding a switch to your choice. top of the receiver; this switch also has an OFF position. When the Caregiver Alert is not in use the receiver should be switched off to extend the battery life.

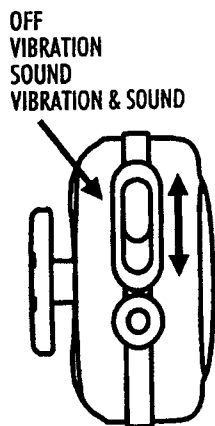

### Function Setting:


Your Caregiver Alert receiver has two alerting methods, one a 6-beep sound and the other is a vibration mode. You can also set the receiver for both functions simultaneously. These functions can be set by sliding the switch from the "Off" position to the mode of your choice.

### Security Code Settings:


Your Caregiver Alert system has up to 128 selectable Security codes available. The code is factory set and does not need to be changed unless you are using more than one system or there is outside interference. Other wireless products may cause interference and the Caregiver Alert may not

operate properly. To change the code, simply follow the instructions below.

1. Open the battery compartment door on both the receiver and transmitter. Remove the batteries.

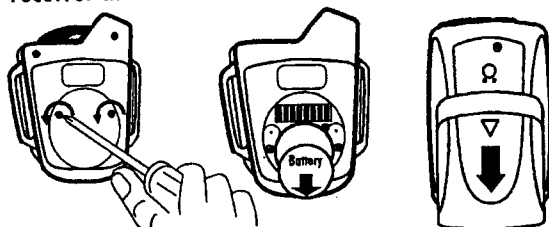

2. Inside both devices you will see a dip switch with switches numbered 1 thru 8.

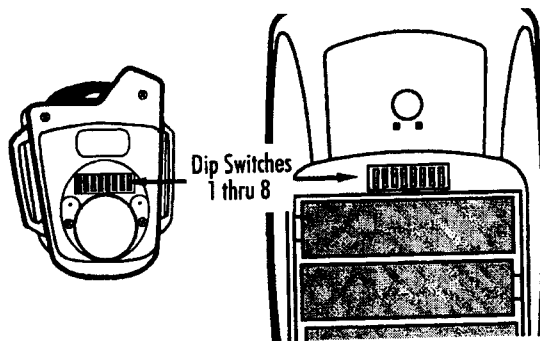

3. Slide any one or more of the switches on the receiver to the ON position. Slide the corresponding switch on the transmitter to the ON position, making sure that both devices have matching settings.

4. Replace the batteries, do not attach the compartment doors. Test system by pressing the button on the transmitter.

5. If system works then replace the battery compartment doors. If it does not recheck the batteries alignment and that the code settings are the same for both units.

### Technical Service

(Do Not Send Products)

If you experience a problem, follow this guide. You may also want to visit our Web site at: [www.desatech.com](http://www.desatech.com). If the problem persists, call\* for assistance at 1-800-858-8501, 7:30 AM to 4:30 PM CST (M-F). You may also write\* to: DESA Specialty Products™  
P.O. Box 90004, 2901 Industrial Drive Bowling Green, Kentucky USA42102  
ATTN: Technical Service Specialty Products

\* If contacting Technical Service, please have the following information available: Model Number, Date of Purchase, and Place of Purchase.

**No Service Parts Available for this Product**

### Regulatory Information

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. The user is cautioned that changes or modifications not expressly approved by the party responsible for regulatory compliance could void the user's authority to operate the equipment.

**Warning: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.**

**NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:**

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

#### ONE YEAR LIMITED WARRANTY

This is a "Limited Warranty" which gives you specific legal rights. You may also have other rights which vary from state to state or province to province. For a period of one year from the date of purchase, any malfunction caused by factory defective parts or workmanship will be corrected at no charge to you. **Batteries are not covered.** To obtain a refund or a replacement, return the product to the place of purchase.

**Not Covered** - Repair service, adjustment and calibration due to misuse, abuse or negligence, light bulbs and other expendable items are not covered by this warranty. Unauthorized service or modification of the product or of any furnished component will void this warranty in its entirety. This warranty does not include reimbursement for inconvenience, installation, setup time, loss of use, or unauthorized service.

This warranty covers only DESA Specialty Products™ assembled products and is not extended to other equipment and components that a customer uses in conjunction with our products.

THIS WARRANTY IS EXPRESSLY IN LIEU OF ALL OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING ANY WARRANTY, REPRESENTATION OR CONDITION OF MERCHANT ABILITY OR THAT THE PRODUCTS ARE FIT FOR ANY PARTICULAR PURPOSE OR USE, AND SPECIFICALLY IN LIEU OF ALL SPECIAL, INDIRECT, INCIDENTAL, OR CONSEQUENTIAL DAMAGES.

REPAIR OR REPLACEMENT SHALL BE THE SOLE REMEDY OF THE CUSTOMER AND THERE SHALL BE NO LIABILITY ON THE PART OF DESA SPECIALTY PRODUCTS™ FOR ANY SPECIAL, INDIRECT, INCIDENTAL, OR CONSEQUENTIAL DAMAGES, INCLUDING BUT NOT LIMITED TO ANY LOSS OF BUSINESS OR PROFITS, WHETHER OR NOT FORESEEABLE. Some states or provinces do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you. Retain receipt for warranty claims.

DESA Specialty Products™ reserves the right to discontinue and to change specifications at any time without notice without incurring any obligation to incorporate new features in previously sold products.

©2005 DESA Specialty, LLC Bronx, NY 10461  
598-1216-01

RECEIVED 11 AUG 2005