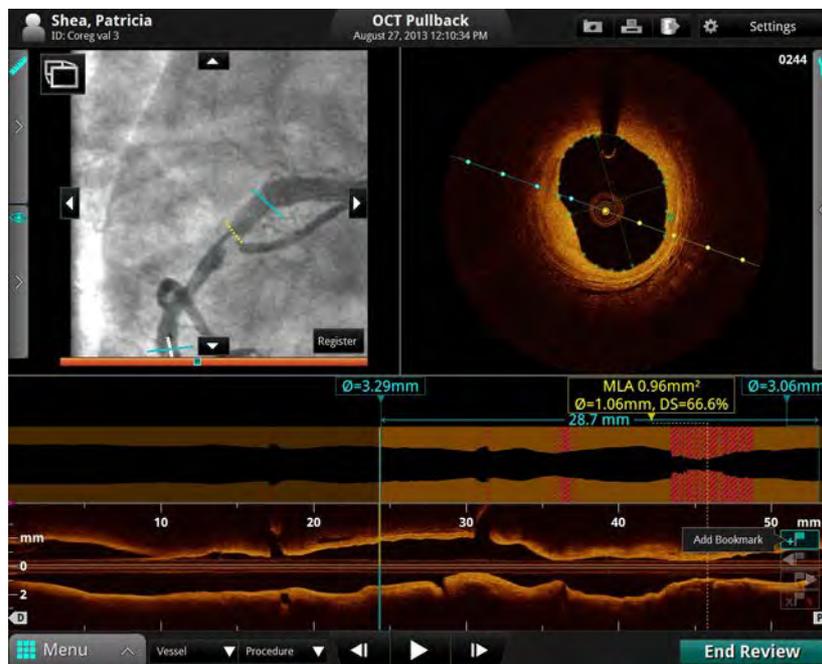


# OPTIS™ Integrated System Instructions for Use



---

Information provided within this Document is subject to change without notice and although believed to be accurate, St. Jude Medical, Inc. and its affiliated companies including without limitation, St. Jude Medical System AB (Sweden), assume no responsibility for any errors, omissions or inaccuracies.

©2014 St. Jude Medical, Inc. and its related companies. All rights reserved. Reproduction, adaptation, or translation without prior permission is prohibited, except as allowed under copyright laws.

OPTIS Integrated systems are subject to US Patents 8,412,312, 8,325,419, 6,565,514; GB, FR 0883793, DE 69738291. ILUMIEN, OPTIS, ST. JUDE MEDICAL, LIGHTLAB IMAGING, GOLDEN IMAGE and the color gold are registered or unregistered trademarks of St. Jude Medical, Inc. and its related companies.

The OPTIS INTEGRATED system software incorporates third party licensed software as described at the following URL: [www.sjmprofessional.com/ilumien-legal-notice](http://www.sjmprofessional.com/ilumien-legal-notice)



St. Jude Medical Coordination Center BVBA  
The Corporate Village  
Da Vincilaan 11 Box F1  
1935 Zaventem Belgium  
+32 2 774 68 11

Australian Sponsor:  
St. Jude Medical Australia Pty Limited  
17 Orion Rd.,  
Lane Cove NSW 2066 Australia  
+61 2 9936 1200



LightLab Imaging, Inc.  
4 Robbins Road  
Westford, MA 01886  
USA

Phone: +1 855 478 5833 US Toll-free  
+1 651 756 5833 International

[www.sjm.com](http://www.sjm.com)

Service E-mail: [OCTservice@sjm.com](mailto:OCTservice@sjm.com)

Part Number ARTUS100109403

ENGLISH

Printed in the U.S.A. 4 / 2014

**CAUTION: Federal law restricts this device to sale by or on the order of a Physician licensed by law of the state in which he practices to use or order the use of the device.**

---

---

## SAFETY INFORMATION

Please review this manual carefully before using your OPTIS Integrated System, especially the safety information in [Chapter 11 “Safety Information”](#). Also, especially note Warnings and Cautions shown throughout the manual.

### WARNINGS

#### **Electrical Shock Hazard**

Do not remove OPTIS INTEGRATED System covers. To avoid electrical shock, use only the power cord supplied with the system and connect only to properly grounded wall outlets. See [Chapter 2 “System Setup”](#) and [Chapter 11 “Safety Information”](#) for electrical safety information.

#### **Explosion Hazard**

Do not operate the OPTIS INTEGRATED System in the presence of flammable anesthetics. Doing so could lead to an explosion.

#### **Visible and Invisible Laser Radiation**

Do not stare into the beam or view the beam directly with optical instruments. Doing so may cause serious eye damage and hazardous radiation exposure.

### **Using These Instructions for Use**

This manual describes the St. Jude Medical OPTIS Integrated System. In it you will find:

- Descriptions of system components and user interface.
- Procedures for using the system for both Optical Coherence Tomography (OCT) and Fractional Flow Reserve Rate (FFR) procedures.
- System safety and cleaning information.
- System specifications.

---

## Conventions Used in These Instructions for Use

- These abbreviations and shortened expressions are used throughout this manual.

Table Front-1: Instructions for Use Conventions

---

Dialog Box or Menu	A window that solicits a user response.
Click or Left-click	Click the left mouse button or the Select button on the Table Side Controller.
Double-click	Click the left mouse button or the Select button on the Table Side Controller twice in quick succession.
Select	Highlight a menu item with the mouse or the tab or arrow keys or the Navigation Controller on the Table Side Controller.

---

- Software text, and physical button labels are shown in bold type. Examples: the **Patient Summary** menu, **Settings** button.
- Keyboard keys are shown enclosed in carets. Examples: <Enter>, <Ctrl>, <Tab>.
- Warnings, Cautions, and Notes are set apart from other text and appear as shown below.

**WARNING:** Warnings alert the user to the possibility of injury, death, or other serious adverse reactions associated with product use or misuse.

**CAUTION:** Cautions alert the user to the possibility of a problem with the product associated with its use or misuse. Problems can include product malfunctions, product failure, and/or damage to the product or damage to other property, or loss of data.

**NOTE:** Notes provide additional information.

## Other Instructions for Use

Details of the imaging catheter are covered in the *Dragonfly Imaging Catheter Instructions for Use* provided with the catheters and are not covered in this manual.

Details of the PressureWire are covered in the *PressureWire Instructions for Use* and are not covered in this manual.

**NOTE:** Additionally, this manual does not provide detailed discussion of the system components, except as they are used with the OPTIS Integrated System.

# Contents

---

## Figures

## Tables

## System Overview

Optis Integrated System Features .....	1-1
OPTIS Integrated System Components .....	1-2
OPTIS Integrated System Accessories .....	1-3
Optis Integrated System Cabinet and Control Room Components .....	1-4
Optis Integrated System - Mobile Workstation .....	1-5
OPTIS Integrated System Table Side Controller .....	1-6
OPTIS Integrated System DOC Holster .....	1-7
OPTIS Integrated Symbols .....	1-8
The Drive-motor and Optical Controller (DOC) .....	1-10
The Wi-Box .....	1-11
Indications for Use and Intended Use .....	1-12
Contraindications .....	1-13
Warnings (OCT) .....	1-14
Precautions (OCT) .....	1-15
Complications (OCT) .....	1-16
Warnings and Precautions (FFR) .....	1-17
Connecting to External Equipment/Accessories .....	1-17
Mechanical Enclosure .....	1-17
Electrical .....	1-17
Electronic Interference .....	1-17

## Contents

---

Aortic Reference Pressure .....	1-18
Pressure Averaging (Mean Setting) .....	1-18
Defibrillation .....	1-18
Recording .....	1-18

## System Setup

Positioning the System .....	2-1
Connecting Your System .....	2-2
System Connections .....	2-2
Powering On and Shutting Down Your System .....	2-2
Power On System Cabinet .....	2-3
Shut Down .....	2-4
FFR Settings .....	2-5
Monitor Setup .....	2-6
Setting Monitor Functions .....	2-6
Setting Monitor Position .....	2-6

## Opening a Patient Record

Select Patient Menu (Home Menu) .....	3-2
Patient Summary Menu .....	3-4
Entering New Patient Information .....	3-6
Editing Patient Information .....	3-8
Editing Case Information .....	3-9
Importing a Patient Database .....	3-10
Opening a Saved Recording or Still Image .....	3-11

## Performing an FFR Procedure

Overview .....	4-1
Required Material and Equipment .....	4-1
Setting up the OPTIS Integrated System .....	4-2
Setting up the Wi-Box with the OPTIS Integrated System .....	4-3
Setting up the PressureWire .....	4-4

---

---

Preparing to Record FFR.....	4-5
Recording FFR.....	4-11
Reviewing an FFR Recording.....	4-13
PressureWire Troubleshooting.....	4-15

## Performing an OCT Procedure

Overview.....	5-1
Required Material and Equipment.....	5-1
OCT Imaging Overview.....	5-2
OCT Operating Modes.....	5-3
OCT Recording Types.....	5-3
OCT Trigger Types.....	5-4
Angio-Coregistration.....	5-5
Setting up the OPTIS Integrated System.....	5-5
Setting up the DOC.....	5-6
Setting up the Dragonfly Imaging Catheter.....	5-7
Preparing to Acquire OCT Recordings.....	5-8
Confirm Recording Settings.....	5-12
Dragonfly Imaging Catheter Insertion and Positioning.....	5-13
Acquiring Patient Images.....	5-16
Removing the Dragonfly Imaging Catheter.....	5-18
Troubleshooting OCT Acquisition.....	5-19
Immediately Stopping DOC Operation.....	5-19
Catheter Failure.....	5-20

## Reviewing OCT Recordings

Image Window.....	6-2
L-Mode View.....	6-4
Limitations of L-Mode Data.....	6-5
Playback Controls.....	6-6
Calibration Adjustment.....	6-7

## Contents

---

Adjust Playback Settings .....	6-9
Bookmark Controls .....	6-10
Setting Playback Range .....	6-11
Exporting a Recording or Still Frame .....	6-12
Capturing Still Images .....	6-12
Saving a Still Image .....	6-12
Printing Still Images .....	6-13
Printing a Still Image .....	6-13
Reviewing with Angio-Coregistration .....	6-14

## Measurements and Annotations

Measurements and Text Callouts in the Image Files .....	7-1
Measurement and Annotation Tools .....	7-2
Verifying Calibration .....	7-3
Techniques to Improve Measurement Accuracy .....	7-3
Measurements and Annotations in the L-Mode View .....	7-3
Length Measurements .....	7-4
Making a Length (Distance) Measurement .....	7-5
Area Measurements .....	7-6
Making a Manual Area Measurement .....	7-6
Adding Text Callouts .....	7-8
Adding Text Callouts .....	7-9
The %AS Calculation .....	7-10
Formula for %AS Calculation .....	7-10
Make a %AS Calculation .....	7-10
The %DS Calculation .....	7-12
Formula for %DS Calculation .....	7-12
Make a %DS Calculation .....	7-12
Field of View .....	7-14
Increase/Decrease Field of View .....	7-14
Zooming In Manually .....	7-15
Quick Zoom .....	7-16
Editing Measurements and Annotations .....	7-17
Moving Individual Points .....	7-18
Adding Points to a Multiple Point Area .....	7-18

---

---

Deleting Points from a Multiple Point Area . . . . .	7-18
Deleting Individual Measurements or Text Callouts . . . . .	7-19
Deleting All Measurements and Text Callouts . . . . .	7-19
Lumen Profile Display Option . . . . .	7-20
Lumen Profile Display With MLA Controls Overview . . . . .	7-21
3D Display Option . . . . .	7-23
3D Tissue Controls . . . . .	7-24
3D Display with Lumen . . . . .	7-25
Limitations of 3D Display . . . . .	7-26

## Exporting, Importing, and Managing Files

Compatible Transfer Media and USB Devices . . . . .	8-2
Optical Media . . . . .	8-2
USB Connected Media . . . . .	8-2
File Formats . . . . .	8-4
About Native (Raw OCT) Format . . . . .	8-4
About DICOM Format . . . . .	8-4
About Standard Format . . . . .	8-5
Image Format and Size in Standard Formats . . . . .	8-6
File Size . . . . .	8-6
Standard File Format . . . . .	8-6
Exporting Files During a Review . . . . .	8-7
Exporting Files in Native (Raw) Format . . . . .	8-8
Exporting Files in DICOM Format . . . . .	8-10
Exporting Files in Standard Formats . . . . .	8-12
Exporting Files from the Patient Summary Menu . . . . .	8-14
Using Exported Standard Format Recordings . . . . .	8-16
St. Jude Medical DICOM Viewer . . . . .	8-16
Importing Files from a CD/DVD or USB . . . . .	8-20
Deleting Files . . . . .	8-22
Deleting Files from the Patient Summary Menu . . . . .	8-22
Deleting Files from the Database Menu . . . . .	8-23
Transfer and Import Messages . . . . .	8-24
Duplicate File Name Messages . . . . .	8-26
Database Statistics . . . . .	8-26
Importing Patient Information From a DICOM Worklist or Storage Server . . . . .	8-27

## Cleaning & Maintenance

Contacting St. Jude Medical Service . . . . .	9-2
Cleaning . . . . .	9-2
Routine Cleaning Procedure . . . . .	9-3
Maintenance . . . . .	9-4
Optical Connection Cleaning Procedure . . . . .	9-5
Optical Adapter Replacement Procedure . . . . .	9-8
Cable Connection Inspection Procedure . . . . .	9-10
Transferring Log Files . . . . .	9-11
Identifying the Software Version . . . . .	9-13
Infection Control . . . . .	9-13
User Troubleshooting . . . . .	9-14
System Disposal . . . . .	9-17

## User Interface Reference

Setup Dialog Box and Submenus . . . . .	10-1
Setup - Acquisition Menu . . . . .	10-2
Setup - Acquisition/Other Menu . . . . .	10-4
Setup - Administration Menu . . . . .	10-5
Setup - Database Menu . . . . .	10-7
Setup - Database/Maintenance Menu . . . . .	10-9
Setup - Database/Physician Menu . . . . .	10-11
Setup - DICOM Menu . . . . .	10-12
Setup - DICOM/Image Options Menu . . . . .	10-17
Setup - DICOM/Local Host Menu . . . . .	10-19
Setup - Display Menu . . . . .	10-21
Setup - Measurements Menu . . . . .	10-23
Setup - Measurements/Labels Menu . . . . .	10-24
Setup - Print Menu . . . . .	10-25
Setup - Service Menu . . . . .	10-26
Setup - Service/System Diagnostics Menu . . . . .	10-29

---

**Safety Information**

Patient Safety . . . . . 11-2  
     General . . . . . 11-2  
     Techniques to Minimize Patient Exposure. . . . . 11-2

Operator Safety . . . . . 11-3  
     Avoiding Operator Light Emission Hazards . . . . . 11-3  
     Repetitive Strain Injury (RSI) . . . . . 11-3

Moving the System . . . . . 11-4

Avoiding Electrical Hazards . . . . . 11-5

Making Proper Electrical Connections . . . . . 11-6

Explosion Hazard. . . . . 11-7

System Imaging Limitations . . . . . 11-7  
     Considerations for Optimal Vessel Imaging . . . . . 11-7  
     Considerations for Optimal Tissue Imaging . . . . . 11-7

Electromagnetic Compatibility . . . . . 11-8

Electromagnetic Interference. . . . . 11-8

Safety Functions Built Into the Ilumien Optis System . . . . . 11-9

**System Specifications**

System - Safety & Regulatory . . . . . 12-1

System - Electrical and Physical . . . . . 12-3

Imaging Specifications . . . . . 12-4

FFR Specifications. . . . . 12-5

Electromagnetic Emissions . . . . . 12-6

Electromagnetic Immunity . . . . . 12-7

Recommended Separation Distances. . . . . 12-11

FCC Statement. . . . . 12-12  
     Essential Performance Defined by Operating Mode . . . . . 12-12

**Index**



# Figures

---

1-1	Optis Integrated System Cabinet and Control Room Components . . . . .	1-4
1-1	Optis Integrated System Mobile Workstation . . . . .	1-5
1-2	Wi-Box in cathlab configuration . . . . .	1-11
2-1	Startup Screen. . . . .	2-3
2-2	Shutdown Menu . . . . .	2-4
3-1	Select Patient Menu . . . . .	3-2
3-2	Default Patient Alert. . . . .	3-3
3-3	Patient Summary Menu . . . . .	3-4
3-4	Add New Patient Menu . . . . .	3-6
3-5	Edit Patient Menu. . . . .	3-8
3-6	Case Information Menu . . . . .	3-9
3-7	Recording as shown in the Patient Summary Menu. . . . .	3-11
3-8	OCT Review Screen. . . . .	3-12
4-1	Cathlab with FFR. . . . .	4-1
4-2	Room Manager: Enter Room Information . . . . .	4-3
4-3	Set AO transducer height and open AO transducer guidance message . . . . .	4-6
4-4	Flush PressureWire guidance message. . . . .	4-7
4-5	Turn on PressureWire guidance message. . . . .	4-8
4-6	Advance PressureWire and Equalize guidance message . . . . .	4-9
4-7	Pd/Pa waveforms equalizing . . . . .	4-11
4-8	Recording . . . . .	4-12
5-1	DOC Connections. . . . .	5-6
5-2	Purge Catheter guidance message . . . . .	5-9
5-3	Plug Catheter into DOC Guidance Message . . . . .	5-10
5-4	Dragonfly Catheter Connected to the DOC . . . . .	5-10
5-5	Catheter Connected, Initial Calibration done. . . . .	5-11
5-6	OCT Settings Menu (during Recording) . . . . .	5-12
5-7	Incorrect and Correct Calibration (Dragonfly Duo shown in fingertips) . . . . .	5-14
5-8	Incorrect and Correct Calibration (Dragonfly OPTIS shown in fingertips). . . . .	5-14
5-9	System Display - Acquisition. . . . .	5-16
5-10	Catheter Failure message . . . . .	5-20
5-11	Safe Unload Guidance, Screen 1 . . . . .	5-21
5-12	Safe Unload Guidance, Screen 2 . . . . .	5-21

## Figures

---

5-13	Safe Unload Guidance, Screen 3	5-22
6-1	Playback Calibration (DragonFly Duo shown, in progress)	6-7
6-2	Field of View Settings	6-9
6-3	Adjusted Playback Range	6-11
6-4	Angio-Coregistration Thumbnail Image	6-14
6-5	View Menu, Angio Co-Registration Button	6-15
6-6	Angio Co-Registration: Register	6-16
6-7	Angio Co-Registration Guidance, Step 1	6-16
6-8	Angio Co-Registration Guidance, Step 2	6-17
6-9	Angio Co-Registration Guidance, Step 2 with Trace	6-18
6-10	Angio Co-Registration Guidance, Step 3	6-19
6-11	Angio Co-registration completed successfully screen	6-20
6-12	Main Screen showing Angio Co-Registration	6-20
7-1	Tools	7-2
7-2	Length Measurement	7-4
7-3	Manual Area Measurement (in progress)	7-6
7-4	Text Callouts	7-8
7-5	Enter Note Dialog Box	7-9
7-6	Select Area Measurement Dialog Box	7-10
7-7	%AS Calculation	7-11
7-8	%AS Error Message	7-11
7-9	Select Length Measurement Dialog Box	7-12
7-10	%DS Calculation	7-13
7-11	Zooming an Image	7-15
7-12	Setup Menu, Display Option	7-16
7-13	3D Display with Lumen	7-25
7-14	MLA Frames in 3D	7-25
8-1	Export Button (OCT)	8-7
8-2	The Export Wizard - Step 1	8-8
8-3	Define Alternate Patient ID Menu	8-8
8-4	Highlighted Records	8-14
8-5	St. Jude Medical DICOM Viewer - Image View	8-16
8-6	St. Jude Medical DICOM Viewer - Attributes View	8-17
8-7	Import Database Menu	8-20
8-8	Deletion Warning Alert	8-23
8-9	Add Patient - Step 1 (Worklist)	8-27
8-10	Add Patient - Step 2	8-28
8-11	Add Patient - Step 1 (Storage Server)	8-30
9-1	Inserting Cleaner Into Optical Adapter	9-6
9-2	Inserting Cleaner Into Catheter	9-7
9-3	Proper Gripping of Adapter for Removal	9-8
9-4	Alignment of Optical Adapter with Optical Carriage	9-9
9-5	Transfer Event Log Files Menu	9-11

---

---

9-6	System Startup Window . . . . .	9-13
10-1	Setup - Acquisition Menu . . . . .	10-2
10-2	Setup - Acquisition/Other Menu . . . . .	10-4
10-3	Setup - Administration Menu . . . . .	10-5
10-4	Setup - Database Menu . . . . .	10-7
10-5	Setup - Database/Maintenance Menu . . . . .	10-9
10-6	Setup - Database/Physician Menu . . . . .	10-11
10-7	Setup - DICOM Menu . . . . .	10-13
10-8	Configure DICOM Menu . . . . .	10-14
10-9	Setup - DICOM/Image Options Menu . . . . .	10-17
10-10	Setup - DICOM/Local Host Menu . . . . .	10-19
10-11	Setup - Display Menu . . . . .	10-21
10-12	Setup - Measurements Menu . . . . .	10-23
10-13	Setup - Measurements/Labels Menu . . . . .	10-24
10-14	Setup - Print Menu . . . . .	10-25
10-15	Setup - Service Menu . . . . .	10-26
10-16	Setup - Service/System Diagnostics Menu . . . . .	10-29
11-1	Connector Panel Laser Safety Labels . . . . .	11-3
11-2	Electrical Label . . . . .	11-6

## Figures

---

# Tables

---

Front-1	Instructions for Use Conventions.....	Front-iv
1-1	Symbols Description .....	1-8
1-2	DOC Controls.....	1-10
3-1	Select Patient Menu functions .....	3-2
3-2	Patient Summary Menu functions .....	3-4
4-1	FFR Review Screen .....	4-13
5-1	System Display Description - Acquisition.....	5-16
6-1	OCT Display Overview .....	6-2
6-2	L-Mode view .....	6-4
6-3	Playback Controls.....	6-6
6-4	Bookmark Controls .....	6-10
7-1	Measurement and Annotation Tool Functions.....	7-2
7-2	MLA Controls .....	7-21
7-3	Lumen Profile %AS and %DS Calculations .....	7-22
7-4	3D Tissue Controls.....	7-24
8-1	Optical Media Characteristics .....	8-2
8-2	DICOM File Attributes .....	8-18
8-3	Import Database Menu Options .....	8-21
8-4	Transfer Messages .....	8-24
8-5	Duplicate File Name Messages .....	8-26
9-1	User Troubleshooting Tips.....	9-14
10-1	Setup Dialog Box Common Options .....	10-1
10-2	Setup - Acquisition Menu Settings.....	10-2
10-3	Setup - Acquisition/Other Menu Settings .....	10-4
10-4	Setup - Administration Menu Settings.....	10-5
10-5	Setup - Database Menu Settings.....	10-7
10-6	Setup - Database/Maintenance Menu Settings.....	10-9
10-7	Setup - Database/Physician Settings .....	10-11
10-8	Setup - DICOM Menu Settings .....	10-13
10-9	Setup - Configure DICOM Menu Settings.....	10-15

## Tables

---

10-10	Setup - DICOM/Image Options Menu Settings . . . . .	10-17
10-11	Setup - DICOM/Local Host Menu Settings . . . . .	10-19
10-12	Setup - Display Menu Settings . . . . .	10-21
10-13	Setup - Measurements Menu Settings . . . . .	10-23
10-14	Setup - Measurements/Labels Menu Settings . . . . .	10-24
10-15	Setup - Print Menu Settings . . . . .	10-25
10-16	Setup - Service Menu Settings . . . . .	10-27
10-17	Setup - Service/System Diagnostics Menu Settings . . . . .	10-29
12-1	System Safety & Regulatory Specifications . . . . .	12-1
12-2	System Electrical and Physical Specifications . . . . .	12-3
12-3	Imaging Specifications . . . . .	12-4
12-4	FFR Specifications . . . . .	12-5
12-5	Guidance and Manufacturer's Declaration - Electromagnetic Emissions . .	12-6
12-6	Guidance and Manufacturer's Declaration - Electromagnetic Immunity . .	12-7
12-7	Guidance and Manufacturer's Declaration - Electromagnetic Immunity . .	12-9
12-8	Recommended separation distances between portable and mobile RF commu- nications equipment and the Ilumien Optis System	12-11

## OPTIS Integrated System Features

Optical Coherence Tomography (OCT) is an imaging modality that uses fiber-optic technology. The OPTIS INTEGRATED System uses optical imaging catheters that emit near-infrared light to produce high-resolution real-time images. The frequency and bandwidth characteristics of the near-infrared light used in these systems result in image resolution that is superior to typical medical ultrasound images.

Fractional Flow Reserve (FFR) is the ratio of distal coronary arterial pressure to aortic pressure, measured during hyperemia. It provides the maximal blood flow in the presence of a stenosis as a fraction of the achievable blood flow that would exist in the hypothetical situation that the stenosis was not present. The physician may use the FFR parameter, along with knowledge of patient history, medical expertise and clinical judgment to determine if therapeutic intervention is indicated. This functionality is achieved when the OPTIS INTEGRATED System is used in conjunction with the manufacturer's wireless distal intracoronary pressure transducer and a proximal aortic pressure transducer.

The OPTIS Integrated System is built into the catheter lab so that OCT and FFR are immediately available without the need to find, connect, position, and power-on a mobile console. The system allows either the sterile operator or non sterile operator to control system functions during image and FFR acquisition and review, and to view the OCT and FFR images on the main catheter lab monitor boom or on the Mobile Workstation monitor. In addition, the system incorporates angio co-registration, which allows the user to visualize the position of OCT image data on angiography images, tightening the linkage between anatomical assessment with OCT and subsequent therapeutic actions.

**CAUTION:** Medical personnel who use the OPTIS INTEGRATED System must be aware of the system's limitations. Only trained operators can determine if use of the OPTIS INTEGRATED System is appropriate. Be sure to read [Chapter 11 "Safety Information"](#), before operating the OPTIS INTEGRATED System for the first time.

## OPTIS INTEGRATED System Components

The OPTIS INTEGRATED System includes the following components, integrated into a catheter lab:

- A system Cabinet which includes an isolation transformer, a laser imaging engine, and a computer.
- A keyboard, and a mouse.<sup>1</sup>
- A Drive-motor and Optical Controller (DOC).
- A Table Side DOC Holster
- A Table Side Controller
- Aortic pressure and PressureWire receivers.
- A Mobile Workstation which includes a monitor, keyboard, and mouse.

**NOTE:** Use only the power cable and accessories provided with the system. Use of other cables or accessories may negatively affect EMC performance.

**NOTE:** FFR procedures require you to have a Wi-Box installed in your cathlab. See your St. Jude Medical service representative for more information.

**CAUTION:** The above components are integral parts of the OPTIS INTEGRATED System. The hardware and software must not be modified in any way by the customer. Making such modifications may interfere with correct operation and will void system warranties. See your St. Jude Medical service representative for more information.

**NOTE:** Contact your St. Jude Medical service representative whenever there is new construction in catheter lab.

See [Chapter 12 “System Specifications”](#) for more information on system components.

---

1. Brands and models of components may vary from those shown in this manual.

## **OPTIS INTEGRATED System Accessories**

- PressureWire
- Wi-Box
- C7 Dragonfly™ imaging catheter, Dragonfly™ Duo imaging catheter, or Dragonfly™ OPTIS™ imaging catheter.
- Sterile DOC Cover
- 3 ml Syringe

**NOTE:** See your sales representative for order numbers of accessories in your market.

**System Overview**  
**OPTIS Integrated System Components**

---

**OPTIS Integrated System Cabinet and Control Room Components**

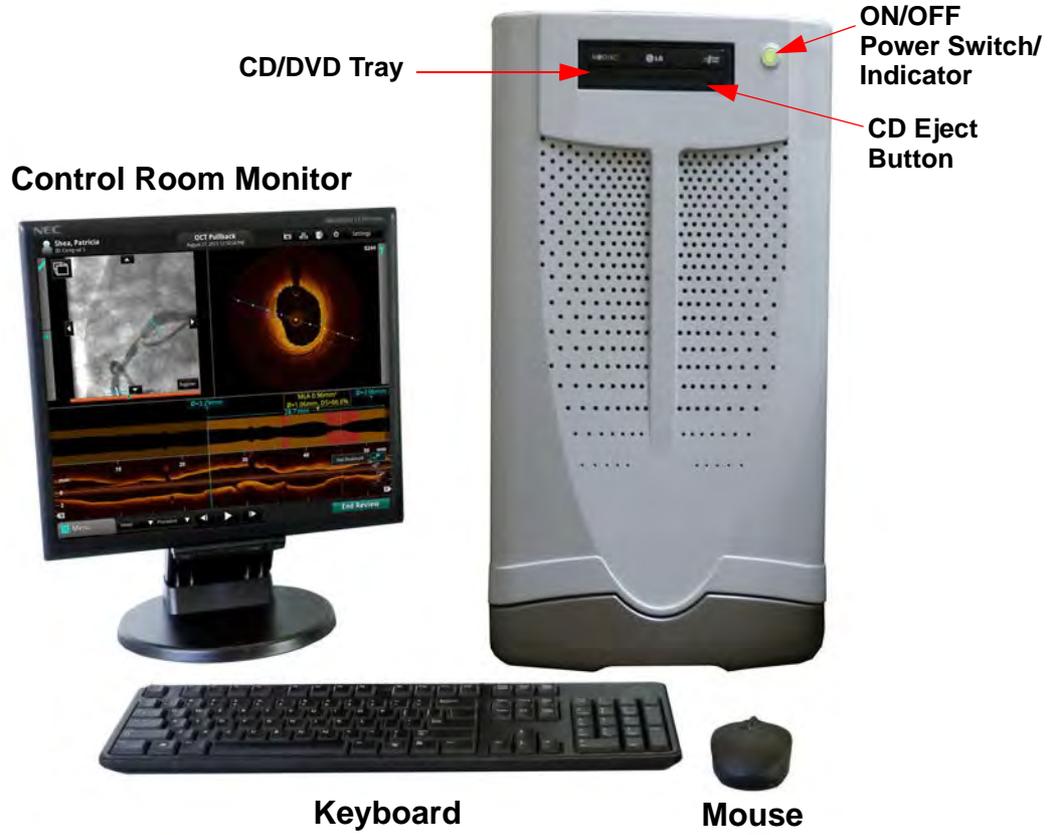


Figure 1-1: OPTIS Integrated System Cabinet and Control Room Components

OPTIS INTEGRATED System - Mobile Workstation

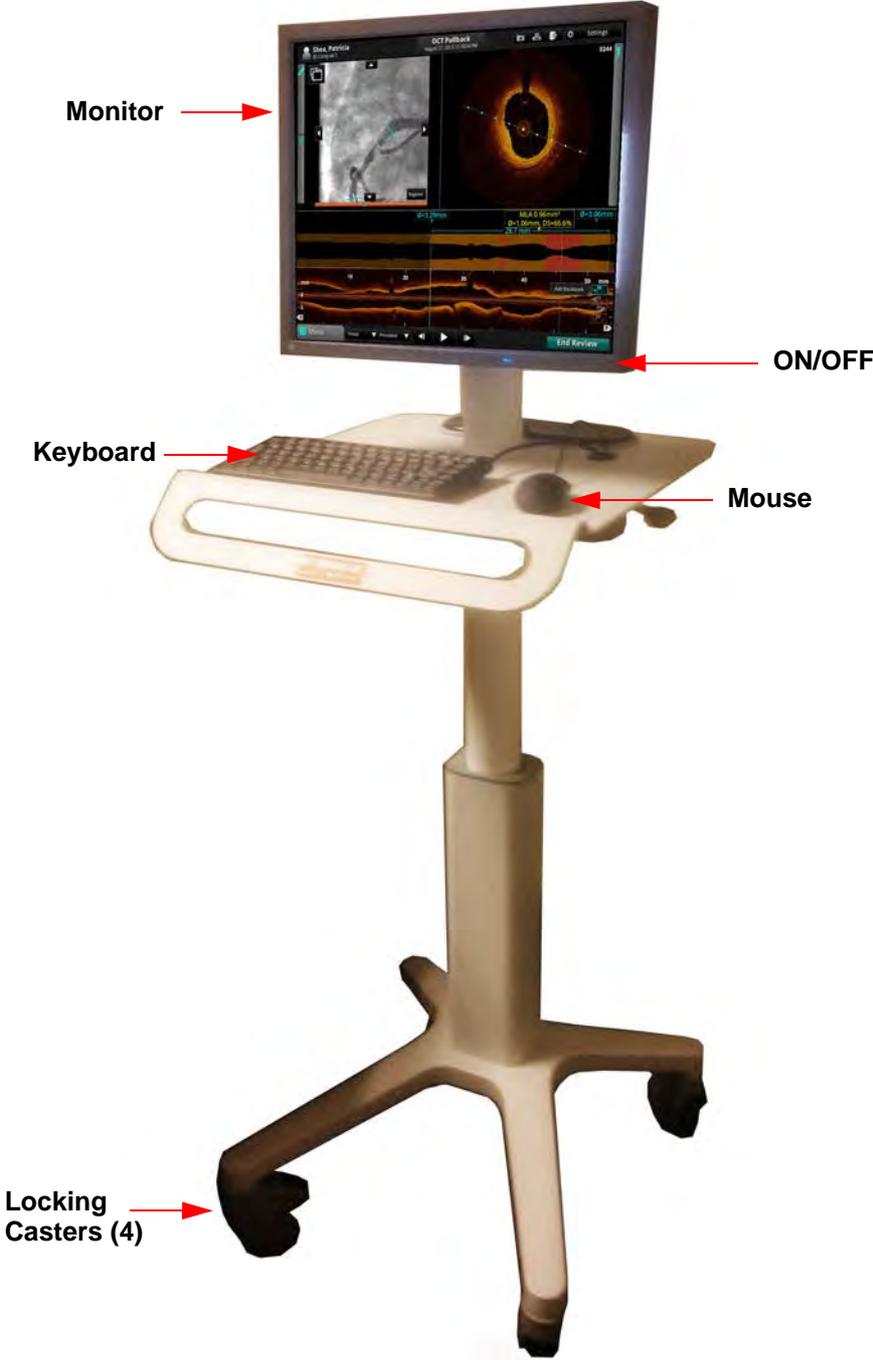


Figure 1-1: OPTIS Integrated System Mobile Workstation

## OPTIS Integrated System Table Side Controller




**Table Side Controller:** The TSC may be used to position the mouse cursor by moving the **Navigation Controller** to left, right, up, or down. The **Select** button (on the top of the Navigation Controller) functions as the left mouse button, and a twist of the **Navigation Controller** functions as the mouse scroll wheel (twist left is scroll up and twist right is scroll down).

**Wireless Connectivity Indicator:**

Blue = Connected

Yellow = No Connection

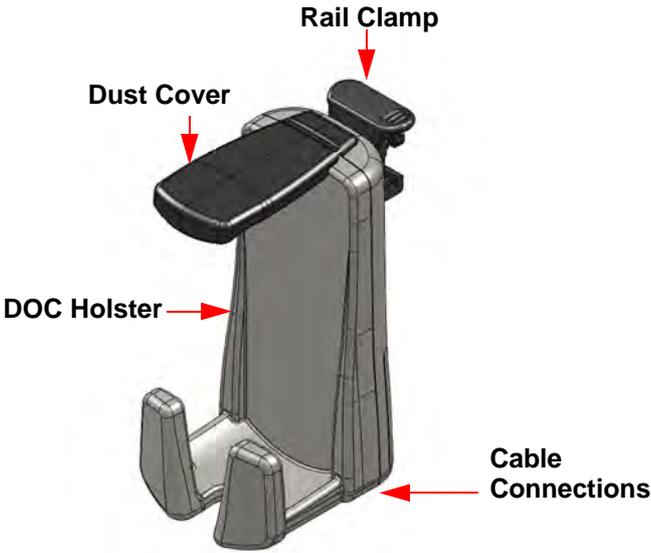**OPTIS Integrated System DOC Holster**



## OPTIS INTEGRATED Symbols

Before using the system, read these Instructions for Use carefully, including the identification of symbols used on the equipment.

Table 1-1: Symbols Description

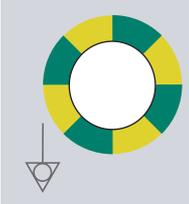
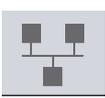
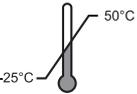
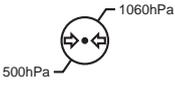
	Equipotential Grounding Post - For a secondary ground connection between equipment.
	European Conformance, BSI Notified Body (British Standards Institution)
	Indoor use only. No protection against ingress of water.
	Curtis Straus NRTL
	Laser hazard symbol - marks a device which produces visible and invisible laser radiation.
	<b>VIDEO:</b> For connection to an external boom monitor.
	<b>ATTENTION!:</b> consult accompanying documents.
	<b>NETWORK:</b> For connecting to a Remote DICOM Server through a network (see the DICOM Setup menus in <a href="#">Chapter 10 “User Interface Reference”</a> ).
	<b>USB:</b> For connection of external storage devices (see <a href="#">“USB Connected Media”</a> on page 8-2).
	<b>WARNING: If the USB device is being used in the patient vicinity it must be port powered or it may compromise electrical isolation and cause patient injury.</b>
	Type CF, defibrillator-proof. Suitable for all patient applications including direct cardiac applications.
	Do not push the system from the monitors or monitor support mount. See <a href="#">“Moving the System”</a> on page 11-4 for all safety instructions on moving the system.

Table 1-1: Symbols Description (*continued*)


	Date of Manufacture
	FCC Identification Number
<del>FCC ID: SB6C408650</del>	
<del>IC: 10934A-101</del>	Industry Canada certification 
 N 17554	Australian Communications and Media Authority supplier code number
	Consult instructions for use
	Interference may occur in the vicinity of this equipment
	Disposal of the equipment must be in accordance with local laws.
	Temperature Range (Shipping label)
	Atmospheric Range (Shipping label)
	Humidity Range (Shipping label)

**WARNING:** All connections to the OPTIS INTEGRATED System must be made through the System Connector Panel. Making connections directly to internal components of the system may bypass isolation features and compromise patient safety.

## The Drive-motor and Optical Controller (DOC)

The Drive-motor and Optical Controller (DOC) provides bed-side control of the most important OCT imaging functions. Refer to “[Preparing to Acquire OCT Recordings](#)” on page 5-8 for Dragonfly Imaging Catheter connection details.

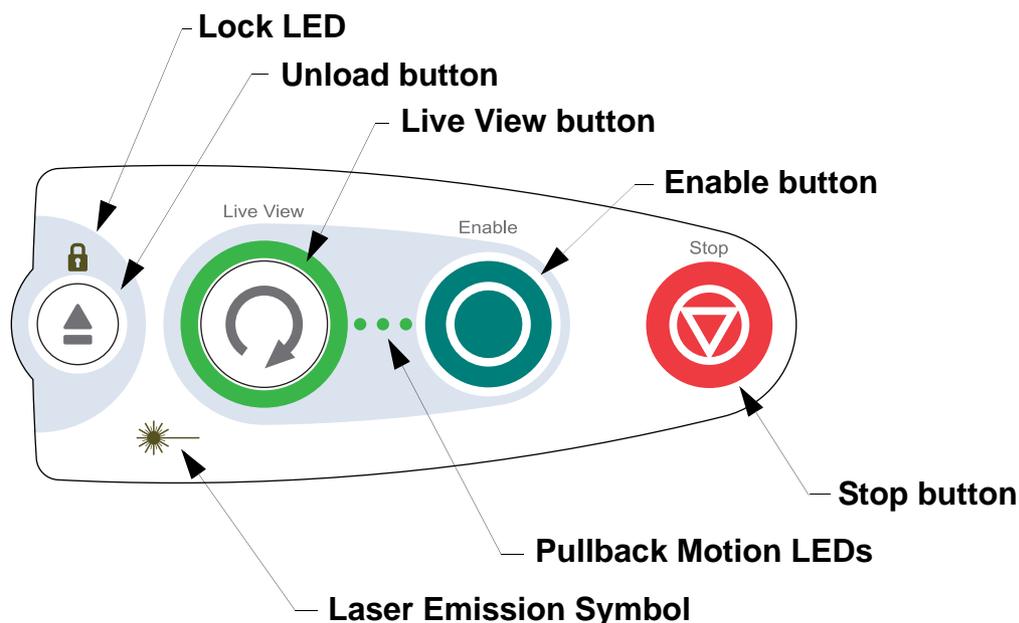Table 1-2: DOC Controls



<b>Lock LED</b>	<ul style="list-style-type: none"> <li>• <b>Off</b> when the imaging catheter is not loaded.</li> <li>• <b>Blinking</b> when loading or unloading the imaging catheter.</li> <li>• <b>On</b> when the imaging catheter is loaded.</li> </ul>
<b>Unload</b>	Press to unload the imaging catheter.
<b>Live View</b>	<p>Press to switch between Standby View and Live View.</p> <p>When the button is lit in green, the system is in Live View.</p> <p>When the button is blinking, the system is enabled to begin a pullback recording.</p> <p><b>NOTE:</b> If the <b>Recording Type</b> is set to <b>Stationary</b>, this button does not blink.</p>
<b>Enable</b>	Press to <b>Auto Calibrate</b> or to <b>Enable</b> recording (if calibration is complete).
<b>Stop</b>	Press to stop imaging catheter motion and turn off laser output.
Pullback Motion LEDs	<ul style="list-style-type: none"> <li>• <b>Off</b> when the imaging catheter is stationary.</li> <li>• <b>Blinks</b> during pullback.</li> </ul>
Laser Emission Symbol	Illuminated when laser output is on (whenever the system is in Live scanning mode).

## The Wi-Box

The Wi-Box is installed in your cathlab between your Hemodynamic Recording System and the AO Transducer. The position of the Wi-Box in a cathlab is shown in Figure 1-2.

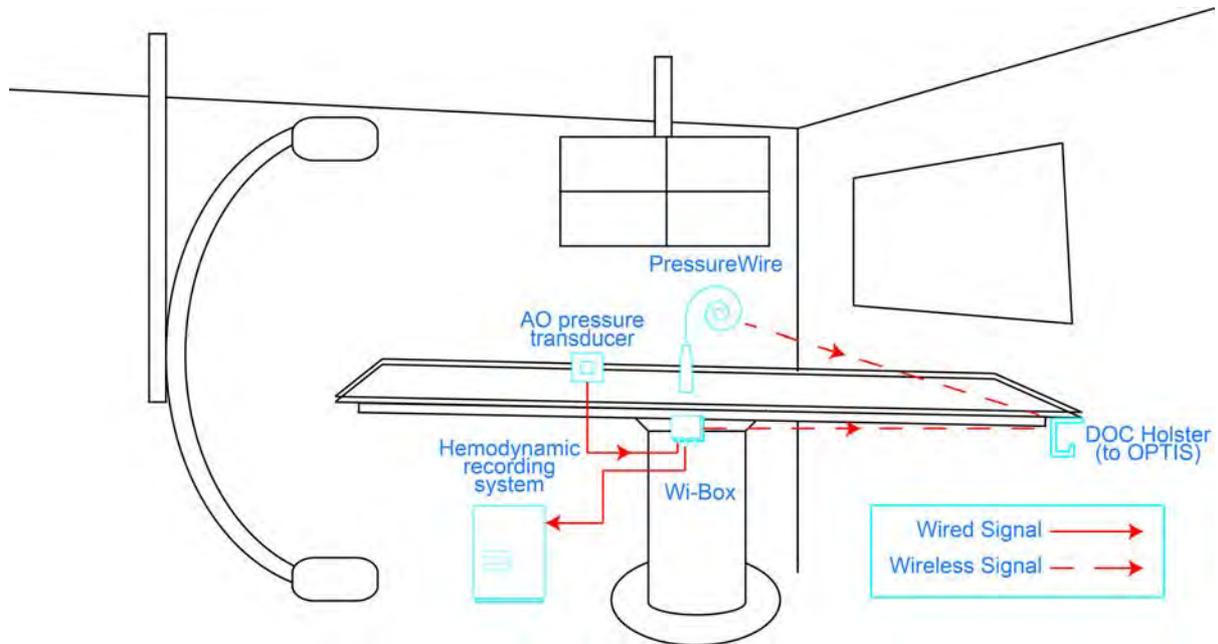


Figure 1-2: Wi-Box in cathlab configuration

The wireless connection to your Wi-Box is made during setup for your procedure. See “Setting up the Wi-Box with the OPTIS Integrated System” on page 4-3 for more information.

## Indications for Use and Intended Use

The OPTIS INTEGRATED Imaging System with Dragonfly Imaging Catheter is intended for the imaging of coronary arteries and is indicated in patients who are candidates for transluminal interventional procedures. The Dragonfly Imaging Catheter is intended for use in vessels 2.0 to 3.5 mm in diameter. The Dragonfly Imaging Catheter is not intended for use in the left main coronary artery or in a target vessel which has undergone a previous bypass procedure.

The OPTIS INTEGRATED Imaging System will further acquire radio frequency signal outputs from both a distal intracoronary pressure transducer and a proximal aortic pressure transducer to determine the physiological parameter, Fractional Flow Reserve (FFR). The physician may use the FFR parameter, along with knowledge of patient history, medical expertise and clinical judgment to determine if therapeutic intervention is indicated.

**CAUTION:** The OPTIS INTEGRATED System is intended for use by appropriate medical personnel who have received OPTIS INTEGRATED training. St. Jude Medical and its employees cannot give instructions in the interpretation or diagnosis of recordings and makes no attempt to do so.

**WARNING:** Prior to use, please review the *Instructions for Use* supplied with the Dragonfly Imaging Catheter and with the PressureWire for more information.

## Contraindications

Use of the St. Jude Medical OPTIS INTEGRATED System is contraindicated where introduction of any catheter would constitute a threat to patient safety.

Contraindications (listed alphabetically) include:

- Bacteremia or sepsis
- Major coagulation system abnormalities
- Patients diagnosed with coronary artery spasm
- Patients disqualified for CABG surgery
- Patients disqualified for PTCA
- Severe hemodynamic instability or shock
- Total occlusion
- Large thrombus
- Acute renal failure

**NOTE:** The system has no patient alarm functions. Do not use for cardiac monitoring.

## **Warnings (OCT)**

- Appropriate anticoagulant and vasodilator therapy must be used during the procedure as needed.
- Observe all advancement and movement of the Dragonfly Imaging Catheter under fluoroscopy. Always advance and withdraw the catheter slowly. Failure to observe device movement fluoroscopically may result in vessel injury or device damage.
- Leave the guidewire engaged with the catheter at all times during use. Do not withdraw or advance the guidewire prior to withdrawing the catheter.
- If resistance is encountered during advancement or withdrawal of the Dragonfly Imaging Catheter, stop manipulation and evaluate under fluoroscopy. If the cause of resistance cannot be determined or mitigated, carefully remove the catheter and guidewire together.
- The catheter should never be forced into lumens that are narrower than the catheter body or forced through a tight or heavily calcified lesion.
- The catheter should not be advanced through abnormally tortuous anatomy.
- When advancing or retracting a catheter with a minirail tip through a stented vessel, the catheter may engage the stent between the junction of the catheter and guidewire, resulting in entrapment of catheter/guidewire, catheter tip separation, and/or stent dislocation.
- Refer to the contrast media's instructions-for-use for general warnings and precautions relating to use of the contrast media.

**CAUTION:** Before creating an OCT recording, review [Chapter 5 “Performing an OCT Procedure”](#) for additional warnings and cautions.

## **Precautions (OCT)**

- Safety and effectiveness have been established for the following patient population: adult patients undergoing non-emergent percutaneous coronary interventions in lesions with reference vessel diameters between 2.0 to 3.5 mm, which were not located in the left main coronary artery or in a target vessel which has undergone previous bypass procedures.
- All operators must be trained prior to using the OPTIS INTEGRATED System and the Dragonfly Imaging Catheter.
- Only 100% contrast media is approved for human use.
- Store the catheter at ambient temperature in a dry location out of direct sunlight.
- Never attempt to attach or detach a catheter to the DOC while the “lock” LED is lit.
- Do not kink, sharply bend, pinch, or crush the catheter at any time.
- The catheter is for single use only. Do not reuse, re-sterilize, or reprocess.
- The catheter is sterilized by ethylene oxide and is intended for one time use only. Non-pyrogenic. Do not use if the package is opened or damaged.
- After use, the catheter may be a potential biohazard. Handle and dispose of in accordance with accepted medical practice and applicable laws and regulations.
- The catheter has no user serviceable parts. Do not attempt to repair or alter any part of the catheter assembly as provided.

## **Complications (OCT)**

The risks involved in vascular imaging include those associated with all catheterization procedures. The following complications (listed alphabetically) may occur as a consequence of intravascular imaging and may necessitate additional medical treatment including surgical intervention.

- Acute myocardial infarction or unstable angina
- Allergic reaction to the contrast media
- Arterial dissection, injury, or perforation
- Cardiac arrhythmias
- Coronary artery spasm
- Death
- Embolism
- Thrombus formation

## Warnings and Precautions (FFR)

**CAUTION:** Before beginning an FFR procedure, review the *PressureWire Instructions for Use*. Details of the PressureWire are covered in the *PressureWire Instructions for Use* and are not covered in this manual.

**CAUTION:** Patients with potential microvascular dysfunction and borderline FFR values should be interpreted with caution, and management strategies should be guided not only by pressure measurement, but also by possibly supplementary clinical risk stratification and other tests.

### Connecting to External Equipment/Accessories

**WARNING:** When used in the patient environment, all equipment connected to the OPTIS INTEGRATED System must meet the requirements for medical isolation according to the IEC 60601 safety standards. Connection of equipment that does not follow relevant IEC standards (e.g. IEC 60601 series for medical electrical equipment) may lead to patient injury or death.

**CAUTION:** No connections to other systems or components are to be made to the OPTIS INTEGRATED System except through the Connector Panel. No connections are to be made through the Connector Panel except as described in this manual.

In addition, all such combinations of systems shall comply with the standard IEC 60601-1-1, Safety requirements for medical electrical systems. Any person who connects external equipment to the OPTIS INTEGRATED system has formed a medical system and is therefore responsible for compliance of the system with the requirements of IEC 60601-1-1. If in doubt contact a qualified technician. Only PressureWire and Wi-Box are intended to be used with the OPTIS INTEGRATED System wireless receivers.

### Mechanical Enclosure

**WARNING:** Do not use the OPTIS INTEGRATED System if it has been dropped or in another way exposed to mechanical or electrical damage or if liquids have penetrated the housing, or the user or patient may be exposed to electrical shock or faulty readings may appear. Contact your supplier for further action.

**CAUTION:** Ensure that all ventilation holes are unblocked or else system overheating and false readings may occur.

### Electrical

**WARNING:** The mains power remains switched on when the system is in STANDBY mode.

Avoid direct or indirect (e.g. via the operator) conductive connection between other electrical equipment and the OPTIS INTEGRATED System. Conductive connection may cause leakage currents to induce ventricular fibrillation. High frequency surgical equipment must not be used on a patient at the same time as PressureWire and the OPTIS INTEGRATED System.

### Electronic Interference

**CAUTION:** Radio transmitting equipment, cellular phones and strong emission sources such as high frequency surgical equipment shall not be used in the close proximity of the OPTIS INTEGRATED System since this could influence the performance of the device.

### Aortic Reference Pressure

**CAUTION:** Check that the monitor cables and AO adapter delivered with the OPTIS INTEGRATED System interface are compatible with the cathlab system to be used. The Aortic Pressure Transducer (AO) should be in accordance with ANSI/AAMI BP22-1994.

Once the lab monitor system has been zeroed, use only the OPTIS INTEGRATED System to calibrate the aortic pressure transducer and PressureWire.

### Pressure Averaging (Mean Setting)

**CAUTION:** Choice of excessively high number of heartbeats may result in insensitive pressure averaging (noticeable when there is a short hyperemic plateau). Choice of excessively low number of heartbeats may result in pressure averaging, which are overly sensitive to arrhythmia and pressure disturbances. An insensitive or overly sensitive mean average of pressure may result in an incorrect FFR value.

### Defibrillation

**CAUTION:** The OPTIS INTEGRATED System is a CF Class I equipment and protected against the effects of a discharge of a defibrillator. PressureWire readings may be affected by defibrillation. Recalibrate PressureWire after defibrillation use.

### Recording

**CAUTION:** The system may place the point of FFR at the wrong location due to abnormal heart beat or artefact in Pa from flushing the guiding catheter. The responsible physician should confirm that the point selected by the system is a valid point of FFR.

**CAUTION:** If the cursor position has been saved, the FFR value is changed accordingly.

**CAUTION:** Before creating an FFR recording, review [Chapter 4 “Performing an FFR Procedure”](#) for additional warnings and cautions.

**CAUTION:** After use, the catheter may be a potential biohazard. Handle and dispose of in accordance with accepted medical practice and applicable laws and regulations.

## Positioning the System

**WARNING:** Failure to position the system components as described may lead to a system tipping hazard or a pinching hazard, causing possible patient or operator injury and damage to the system.

- Confirm that the OPTIS Integrated System Cabinet is on a level surface inside the Control Room or X-Ray Room. Confirm that the Control Room Monitor, keyboard, and mouse are located on the Control Room Table.

**NOTE:** The USB ports are located on the keyboard.

- Mount the Table Side Controller onto the operating table rail in a location that is comfortable for the physician during the procedure.
- Position the optional OPTIS Integrated Mobile Workstation (if used) at the foot of the patient table with the monitor facing the attending physician.

Care must be taken to ensure the Mobile Workstation is clear of any moving equipment, including the angiography system. It is the responsibility of the attending physician to ensure that collisions do not occur.

**NOTE:** Whenever the Mobile Workstation is used near moving equipment, it is recommended that the wheels remain unlocked to allow the system to roll if it is bumped.

- Ensure the power cord and any other connections to the Mobile Workstation are routed to prevent a tripping hazard. Ensure that the main power switch and power plug can be accessed at any time during the procedure.

## Connecting Your System

### System Connections

To connect the Mobile Workstation to power, plug the power cord into a grounded electrical outlet. There is no Main Power Switch on the Mobile Workstation. The System Cabinet and Monitors are normally always connected to power.

### Powering On and Shutting Down Your System

Make sure the power cords are connected to the System Cabinet, optional Mobile Workstation (if used), and control room monitor and they are plugged into grounded electrical outlets. For detailed information on electrical requirements, see [“Making Proper Electrical Connections”](#) on page 11-6.

**WARNING:** Never use a converter adapter to plug the three-pronged AC plug into a two-pronged, ungrounded wall outlet. Doing so may result in electric shock to the patient or operator and damage to equipment.

## Power On System Cabinet

To power on your System Cabinet:

1. Confirm that all monitors are powered by observing that the appropriate power indicators are lit. See [“Monitor Setup” on page 2-6](#) for details.
2. Press the main power switch located at the upper-right side of the front of the cabinet to turn on system power (see [Figure 1-1 on page 1-4](#)).

The system's Startup screen appears ([Figure 2-1](#)).



Figure 2-1: Startup Screen

**NOTE:** The first time the software runs, the End User Licence Agreement (EULA) displays. You must first read and agree to the EULA by checking the “I agree to the terms of the EULA” checkbox in order to proceed.

**NOTE:** Once the system is turned on and the System Display is being shown, you may need to adjust the monitor’s brightness and contrast as described in [“Monitor Setup” on page 2-6](#).

## Shut Down

**CAUTION:** Do not unplug from AC power or turn off main power until the shutdown is complete, the screens turn black, and the monitor LEDs turn amber. Disconnecting from AC power before the shutdown is complete may damage the system.

**NOTE:** Use the Shutdown button at the top of the screen to shut down the system. Use of the ON/OFF switch is not recommended to shut down the system.

To shut down the system:

1. Click the **Shutdown** button located at the top of the screen. The system's **Shutdown** menu appears (Figure 2-2).

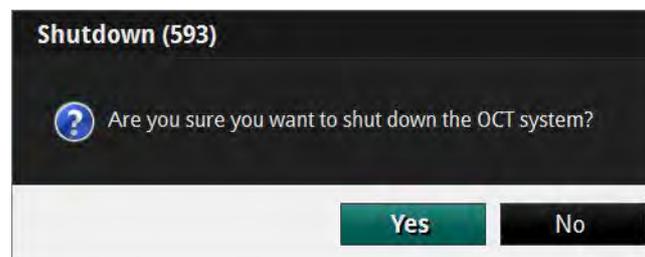


Figure 2-2: Shutdown Menu

2. Click **Yes** to begin the system shutdown, or **No** to continue using the system.

If you choose **Yes**, the computer begins the system shutdown. After 15 seconds, the screens turn black, the monitor LEDs turn amber, and the system enters standby mode.

3. After the screens turn black and the monitor LEDs turn amber, press the main power switch at the base of the cart to turn off system power.
4. If necessary, disconnect the power cord from AC power.

## FFR Settings

**NOTE:** To access the full **FFR Settings** menu, you must be in the process of making an FFR recording, and have selected a room for the FFR procedure.

1. Click the **Settings** button at the top of the screen.

The **FFR Settings** menu opens.

2. Confirm that the FFR settings are correct for this patient.
  - **Pressure Scale:** Changes the vertical scale of the pressure waveform display. The default setting is 0-200 mmHg.
  - **Sweep Speed:** Changes how fast the screen is updated and the level of detail visible to the user. A high number is suitable when a detailed picture of the tracings is required. A low number is suitable when displaying slow changes, for instance during intravenous infusing and pullback. The default setting is **Normal**.
  - **Mean Filter Length (Beats):** Changes the time over which the mean pressure value is calculated. The adjustments are made by selecting the mean calculation filter length, measured in number of heartbeats. The default setting is 3 heartbeats.

**CAUTION:** Choosing a high number of heartbeats makes the pressure averaging slower and less sensitive to artifacts, but may also result in overly insensitive averaging which is noticeable when there is a short hyperemic plateau. Choosing a low number of heartbeats makes the pressure averaging faster and more sensitive to pressure changes, desirable using a short hyperemic plateau, but it may also result in an averaging overly sensitive to arrhythmia and pressure disturbances.

**NOTE:** An insensitive or overly sensitive averaging of pressure may result in an incorrect FFR value.  
Changes in the **Mean Filter Length** setting are not applied to previous FFR recordings.  
When a new patient is selected, the **Mean Filter Length** setting is reset to 3 heartbeats.

## **Monitor Setup**

The flat panel display monitors provided with the System have controls for brightness, contrast, and other monitor functions. Each monitor's display angle is adjustable.

### **Setting Monitor Functions**

Refer to the Instructions for Use that accompanied your monitor.

### **Setting Monitor Position**

The angle of the flat panel monitor should be set to eliminate or minimize glare from surrounding lighting. To optimize viewing:

1. Grasp the monitor firmly with one hand on each side of the screen.
2. Adjust the viewing angle by tilting the monitor.



# Opening a Patient Record

---

3

**CAUTION:** Please note St. Jude Medical makes no representation or warranty that use of the OPTIS INTEGRATED System complies with applicable privacy, security and confidentiality laws, but encourages you to assess your own risk as you use, disclose, control, process or transfer patient health information with the OPTIS INTEGRATED System.

## Opening a Patient Record

### Select Patient Menu (Home Menu)

## Select Patient Menu (Home Menu)

When the system is first started, the **Select Patient** menu is displayed (Figure 3-1). From this menu, you can enter a new patient in the database, open an existing patient, or import a previous patient record.

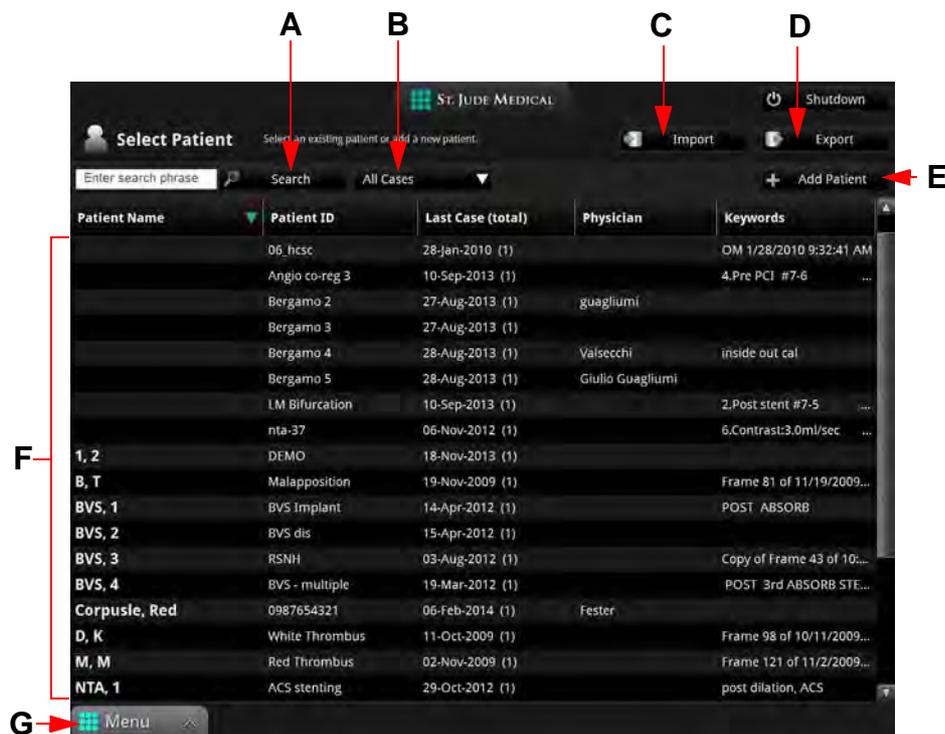


Figure 3-1: Select Patient Menu

Table 3-1: Select Patient Menu functions

<b>A</b>	<b>Search</b> button: Enter a term and click this button to search the patient database.
<b>B</b>	<b>Case Filter</b> drop-down list box: Select <b>ALL</b> cases, <b>OCT</b> only, or <b>FFR</b> only.
<b>C</b>	<b>Import</b> button: Click this button to open the <b>Import</b> menu.
<b>D</b>	<b>Export</b> button: Click this button to open the <b>Export Wizard</b> .
<b>E</b>	<b>Add New Patient</b> button: Click this button to enter a new patient into the System database.
<b>F</b>	List of all existing patients in the system database. Click on a column header to sort the patients according to data in that column. Scroll names with mouse wheel (or twist <b>Navigation Controller</b> ) as desired; single-click name to select patient.
<b>G</b>	<b>Menu</b> : Displays the context-sensitive menu.

**WARNING:** If you select the default patient (“Patient, Default” in the patient list), the system displays an alert (**Figure 3-2**). Do not use the default patient to store patient images. Click OK to continue to use the default patient, or Cancel to return to the Select Patient menu to select another patient.

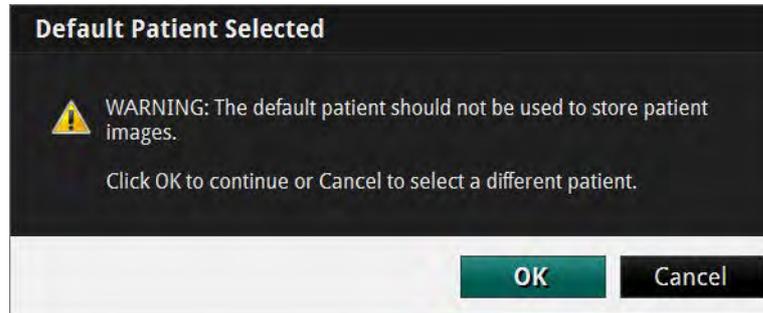


Figure 3-2: Default Patient Alert

## Opening a Patient Record

### Patient Summary Menu

## Patient Summary Menu

When you click on a patient name, the **Patient Summary** menu for that patient opens (see [Figure 3-3](#)). If there are previous recordings for this patient, they are shown here, sorted by date, with the most recent recordings at the top.

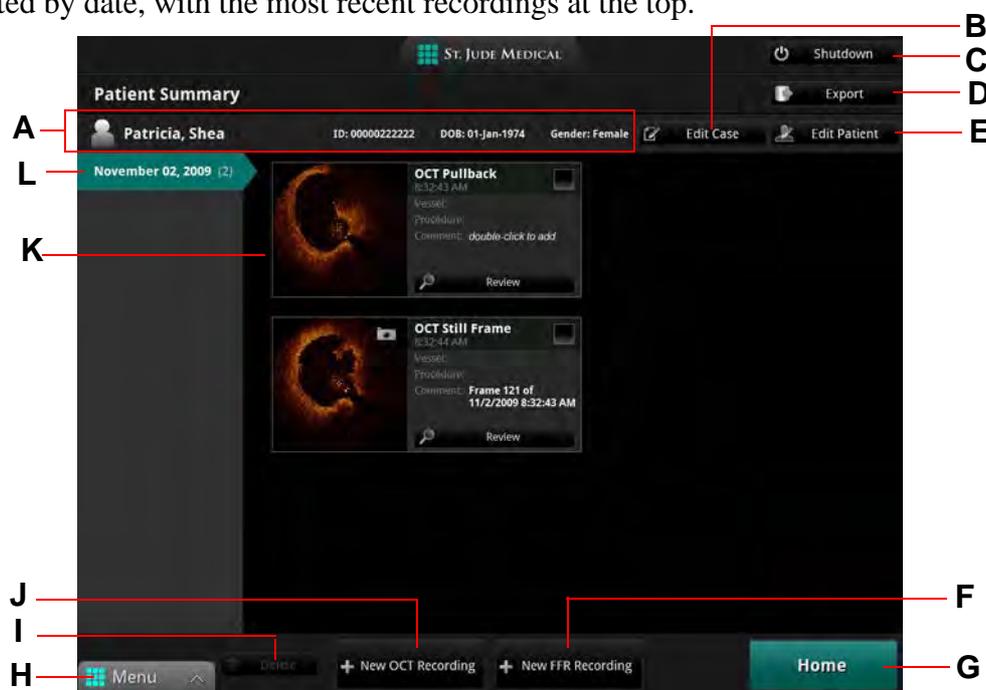


Figure 3-3: Patient Summary Menu

Table 3-2: Patient Summary Menu functions

- |          |   |
|----------|---|
| <b>A</b> | Patient information, including name, ID, date of birth, and gender.   |
| <b>B</b> | <b>Edit Case</b> : Click this button to edit case information.  |
| <b>C</b> | <b>Shutdown</b> : Click this button to begin the shutdown sequence for the system.  |
| <b>D</b> | <b>Export</b> : Click the <b>Checkbox</b> in one or more recordings to select them, then click the <b>Export</b> button to open the <b>Export Wizard</b> .  |
| <b>E</b> | <b>Edit Patient</b> : Click this button to edit patient information.  |
| <b>F</b> | <b>New FFR Recording</b> : Use this button to begin a new FFR recording for this patient.   |
| <b>G</b> | <b>Home</b> : Use this button to return to the <b>Select Patient</b> menu.  |
| <b>H</b> | <b>Menu</b> : Displays the context-sensitive menu. Provides access to the <b>Setup</b> menu.  |
| <b>I</b> | <b>Delete</b> : Check one or more recordings to select them, and click the <b>Delete</b> button to delete them.   |
| <b>J</b> | <b>New OCT Recording</b> : Use this button to begin a new OCT recording for this patient.   |
| <b>K</b> | Gallery of existing recordings for this patient on the date highlighted at left. Double-click in the <b>Comment</b> field to add a comment for that recording. Click the <b>Review</b> button under a recording to open it. |

Table 3-2: Patient Summary Menu functions

---

- L** Case List (dates of recordings) for this patient. The list is sorted by date, with the most recent recordings at the top.
-


## Entering New Patient Information

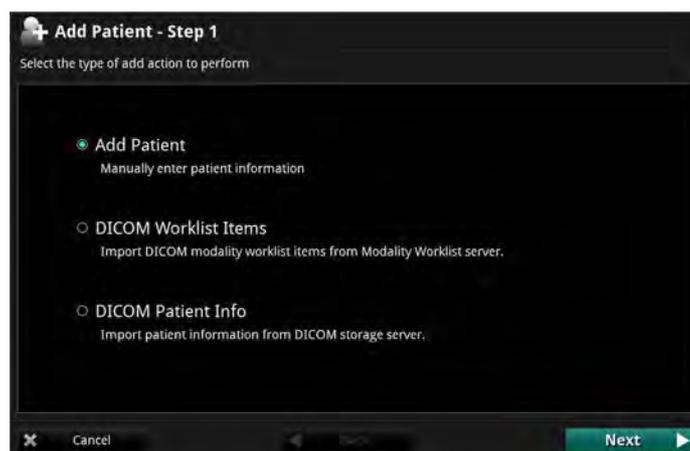You can enter a new patient from the **Select Patient** menu.

1. In the **Select Patient** menu, select the **Add Patient** button. Click **Next**.

The **Add Patient** Wizard displays.

**NOTE:** You can also import patient information from a DICOM storage server or worklist. See [“Importing Patient Information From a DICOM Worklist or Storage Server”](#) on page 8-27.



2. Select **Add Patient**, then click **Next**.

The **Add new patient** menu opens (see [Figure 3-4](#)).



Figure 3-4: Add New Patient Menu

3. Enter the patient information as needed.

**NOTE:** You must enter **Patient ID**, **First name**, and **Last name** before you can save the patient information.

4. Click **New OCT Recording** to save and begin a new OCT recording for this patient, click **New FFR Recording** to save and begin a new FFR recording for this patient, or **Cancel** to close the menu without saving and return to the **Select Patient** menu. Click **Back** to return to the **Add Patient** Wizard.

## Editing Patient Information

The patient information can be edited from the **Patient Summary** menu.

1. In the **Select Patient** menu, click on the patient.

The **Patient Summary** menu for that patient opens (see [Figure 3-3](#)).

2. Click the **Edit Patient** button.

The **Edit Patient** menu opens (see [Figure 3-5](#)).



The screenshot shows a dark-themed dialog box titled "Edit Patient" with a pencil icon. Below the title is the instruction "Enter patient details." The form contains the following fields: "Patient ID:" with the value "F0730"; "Last name:" with the value "Shea"; "First name:" with the value "Patricia"; "DOB:" with a date picker set to "21 Jun 1944"; and "Gender:" with radio buttons for "Male" and "Female", where "Female" is selected. At the bottom right, there are "Cancel" and "OK" buttons.

Figure 3-5: Edit Patient Menu

3. Edit the patient information as needed.

**NOTE:** You must enter **Patient ID**, **First name**, and **Last name** before you can save the patient information.

4. Click **OK** to save the changes, or **Cancel** to close the menu without saving.

## Editing Case Information

The physician name and accession number for a case can be edited from the **Patient Summary** menu.

1. In the **Select Patient** menu, click on the patient.

The **Patient Summary** menu for that patient opens (see [Figure 3-3](#)).

2. Click on a case to select it.

**NOTE:** In the OPTIS INTEGRATED System, all recordings and still images from the same date use the same physician and accession number.

3. Click the **Edit Case** button.

The **Case Information** menu opens (see [Figure 3-6](#)).

The screenshot shows a dark-themed dialog box titled "Case Information". Below the title is the instruction "Enter case details.". There are two input fields: "Physician name:" with the text "Smith" and "Accession number:" with the text "555888". A small green note below the accession number field reads "Optional (for STICOM Support)". At the bottom of the dialog are two buttons: "Cancel" with a close icon and "OK" with a checkmark icon.

Figure 3-6: Case Information Menu

4. Edit the **Physician name** and **Accession Number** as needed.
5. Click **OK** to save the changes, or **Cancel** to close the menu without saving.

## Importing a Patient Database

Previous C7 XR, LUMIEN, and OPTIS INTEGRATED OCT and FFR recordings can be imported into the system using the **Import** button on the **Select Patient** menu. For more information on importing patient files or information, see [“Importing Files from a CD/DVD or USB”](#) on page 8-20.

## Opening a Saved Recording or Still Image

Each recording or still image has an entry in the **Patient Summary** menu. The entry includes a timestamp and a thumbnail, and can display the **Vessel**, **Procedure**, and a **Comment**, if added. The thumbnail adds icons to represent the status and content of the file:

-  • A “camera” icon indicates a still frame from an OCT recording.
-  • A “safe” icon indicates that the file has been archived.

To review a saved recording or still image:

1. Open the **Patient Summary** menu:
  - From the **Select Patient** menu, click on a patient’s name to select.
  - From an OCT or FFR recording, click on the **End Review** button at the bottom of the screen.

The **Patient Summary** menu for that patient opens ([Figure 3-3 on page 3-4](#)).

2. Click the date of the recording from the dates at the left of the menu.

All recordings and still images for that patient on that date are displayed.

3. Click on the **Review** button beneath the recording you want to open.

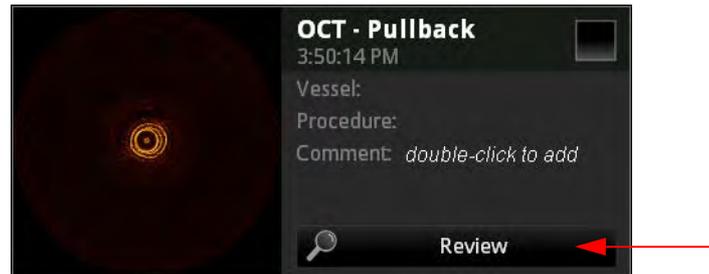


Figure 3-7: Recording as shown in the Patient Summary Menu

The recording or still image opens. Recordings play automatically.

## Opening a Patient Record

### Opening a Saved Recording or Still Image

---



Figure 3-8: OCT Review Screen

- To navigate an OCT recording, see [“Playback Controls”](#) on page 6-6; to navigate an FFR recording, see [“Reviewing an FFR Recording”](#) on page 4-13.
  - To add measurements and annotations to the OCT recording or still frame, see [Chapter 7 “Measurements and Annotations”](#).
  - To export the entire OCT recording or individual frames, see [Chapter 8 “Exporting, Importing, and Managing Files”](#).
4. To end the review and return to the **Patient Summary** menu, click the **End Review** button.

# Performing an FFR Procedure

## 4

### Overview

The FFR procedure requires up to two operators; a sterile operator and optionally, a non-sterile operator. All steps requiring contact with the PressureWire must be performed by the sterile operator. Any steps performed in direct contact with a keyboard or mouse must be performed by the non-sterile operator.

A typical cathlab configured for FFR is shown below.

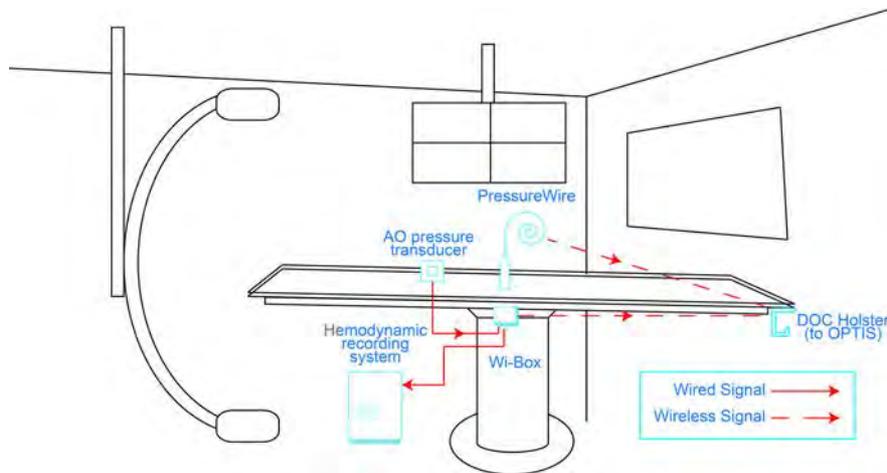


Figure 4-1: Cathlab with FFR

### Required Material and Equipment

- OPTIS INTEGRATED System
- PressureWire
- Wi-Box, installed in each room where the OPTIS INTEGRATED System is to be used to measure FFR
- Heparinized, physiologic saline solution, for hydrophilic catheter preparation

## **Setting up the OPTIS INTEGRATED System**

1. Position the system for use. See [“Positioning the System”](#) on page 2-1.
2. Turn on the system. See [“Power On System Cabinet”](#) on page 2-3.

## Setting up the Wi-Box with the OPTIS INTEGRATED System

The Wi-Box should be connected to your facility's Hemodynamic Recording System at installation. Complete the following steps to connect the OPTIS INTEGRATED System to a Wi-Box for the first time .

1. Verify that the Wi-Box is powered-up by checking the power indicator light. Refer to the installation instructions that came with your Wi-Box for further information.
2. From the Patient Summary screen or from the Select Patient menu, click on the **Menu** button.
3. Click the **Setup** menu option.
4. Click on **Room Manager**. The **Room Manager** displays.

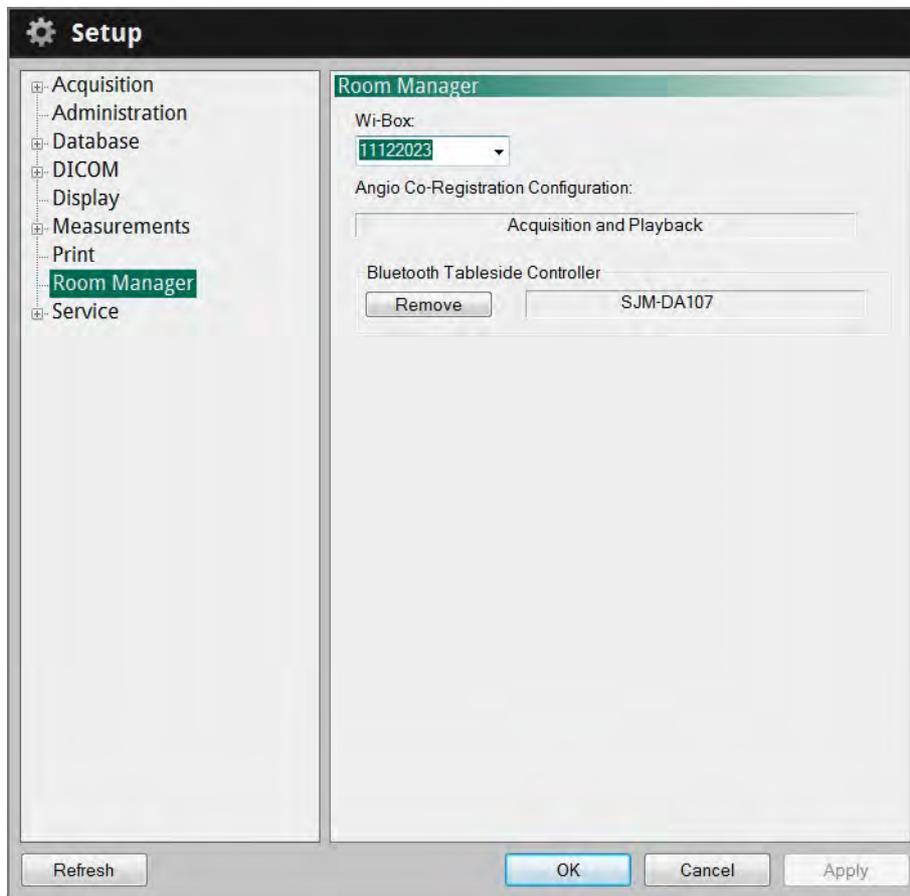


Figure 4-2: Room Manager: Enter Room Information

5. Verify that the Wi-Box Serial Number that appears in this dialog box matches the Serial Number of the Wi-Box in your catheter lab.

## Performing an FFR Procedure Setting up the PressureWire

---

**NOTE:** Be sure to check the Wi-Box serial number. The **Devices** box will auto-detect any active Wi-Box that is in-range. Ensure that the Wi-Box serial number selected matches the Wi-Box device present in the current room. (The serial number is visible on the front panel of the Wi-Box.)

**NOTE:** The Wi-Box serial numbers can also be entered directly into the **Devices** box.

**NOTE:** To change the Wi-Box selection, select the Wi-Box Serial Number in the **Room Manager** dialog box, which enables the **Remove** button. Then, click the **Remove** button to delete the current selection. Then, follow the steps listed in [“Setting up the Wi-Box with the OPTIS Integrated System”](#) on page 4-3 above.

## Setting up the PressureWire

**WARNING:** Please review the *Instructions for Use* supplied with the PressureWire for a complete list of warnings, cautions, and setup instructions.

## Preparing to Record FFR

**WARNING:** The heart rate and mean pressure values shown on the OPTIS INTEGRATED System are for reference only and are not intended to be used as the primary display.

Start FFR recordings from the **Patient Summary** menu.

1. If necessary, enter the patient's information. See [“Entering New Patient Information” on page 3-6](#).

**NOTE:** Patient information should be entered into the system and selected for use before beginning a recording.  
If you are entering a patient through the **Add New Patient** menu, click the **New FFR Recording** button at the bottom of the menu, and continue with [Step 4](#) in this section.

2. In the **Select Patient** menu, click on a patient name.

The **Patient Summary** menu for that patient opens ([Figure 3-3 on page 3-4](#)).

3. In the **Patient Summary** menu, click the **New FFR Recording** button.

**CAUTION:** In order to connect to the correct AO source (Wi-Box) you must select the room where the system is being used.  
The first time you connect to a room, you must enter the room's information into the system. See [“Setting up the Wi-Box with the OPTIS Integrated System” on page 4-3](#) for more information.

## Performing an FFR Procedure Preparing to Record FFR

---

The screen displays the main screen with the guidance message “Set AO transducer height to heart level, then open AO transducer to air. Click Zero Pa.”



Figure 4-3: Set AO transducer height and open AO transducer guidance message

4. Position the AO transducer so that it is level with the patient’s heart.

**NOTE:** The AO transducer should remain level with the patient’s heart throughout the procedure.

5. Open the AO transducer to air.

6. In the Pa measurement box, click the **Zero** button, or on the Table Side Controller press the **Proximal Marker** button.

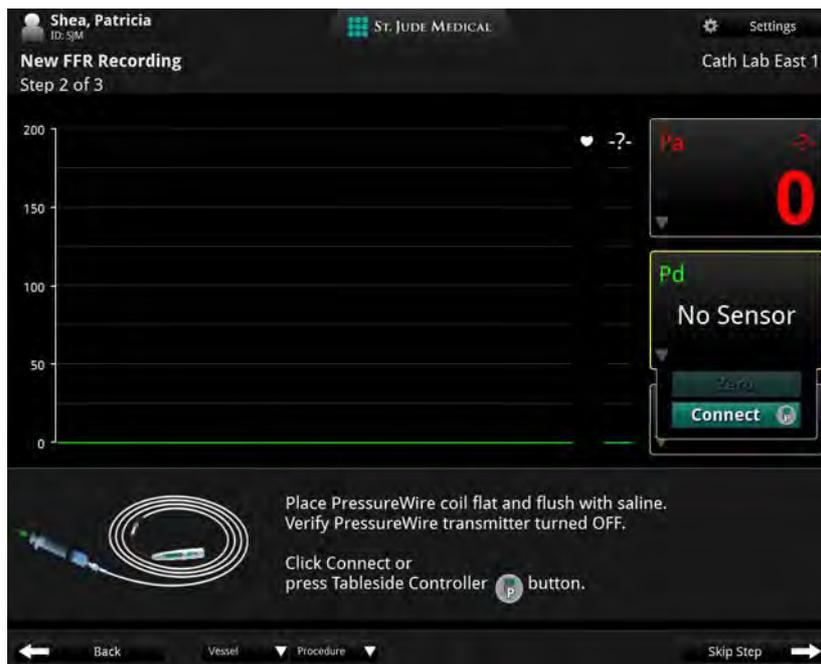


Figure 4-4: Flush PressureWire guidance message

7. Close the AO transducer.
8. Prepare the PressureWire in accordance with the *PressureWire Instructions for Use* and the on screen prompts.

**CAUTION:** Do not use the PressureWire if there are any signs of damage.

9. In the Pd measurement box, click the **Connect** button, or on the Table Side Controller press the **Proximal Marker** button.

**NOTE:** You have 60 seconds to make the connection between the PressureWire transmitter and the system. If necessary, click the **Connect** button again after the 60 second time-out.

## Performing an FFR Procedure Preparing to Record FFR

---

The system is now looking for a PressureWire. The Pd measurement box shows the message “Searching.”

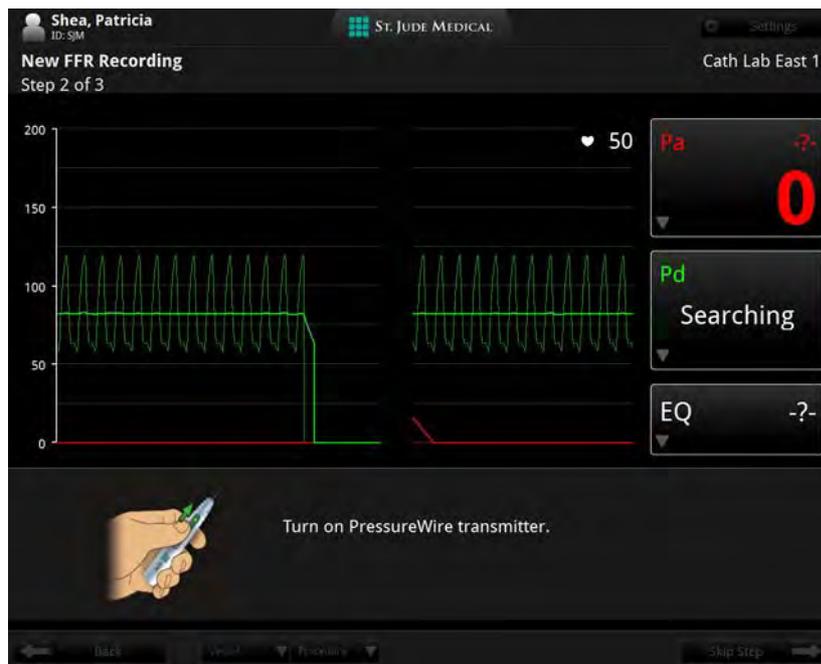


Figure 4-5: Turn on PressureWire guidance message

**CAUTION:** Do not turn on more than one PressureWire while the system is Searching/Connecting.

10. Turn on the PressureWire transmitter.

The Pd measurement box shows the message “Connecting” as the PressureWire transmitter and the system make the wireless connection. Once the connection is established, the Pd waveform appears (in green) in the graphical area of the display.

**NOTE:** If the system does not find the PressureWire after one minute, the message in the Pd measurement box changes back to “No Sensor.” Turn off the PressureWire and return to [Step 9](#).

11. Remove the PressureWire from the plastic hoop.
12. Insert the PressureWire into the patient in accordance with the *PressureWire Instructions for Use*.

The pressure from the AO transducer and PressureWire must now be equalized.



Figure 4-6: Advance PressureWire and Equalize guidance message

**CAUTION:** If the PressureWire is turned off or loses power at any point after equalization, you must remove it from the patient and re-zero it in its plastic hoop in heparinized saline before continuing. See the PressureWire *Instructions for Use* for more information on troubleshooting the PressureWire.

## Performing an FFR Procedure

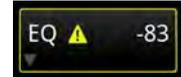
### Preparing to Record FFR

---

13. Check the pressures displayed in the Pa and Pd measurement boxes (See [Figure 4-6 on page 4-9](#)). If the pressures are not equal, click the **Equalize** button (or on the Table Side Controller press the **Proximal Marker** button) to equalize the pressure reading between the AO transducer and the PressureWire.

The equalization offset value appears in the Equalization measurement box. The offset is applied to the Pd pressure, to match the value of the Pa pressure.

**CAUTION:** If the equalization offset value is greater than or equal to 30, or equal to or lower than -30, the system displays an alert icon in the Equalization measurement box. If the system displays this alert icon, try the following to reduce the EQ value:



- Confirm the AO transducer is positioned at the same height as the patient's heart and re-equalize.
- If the alert remains, remove the PressureWire from the patient and re-zero it in its plastic hoop in heparinized saline.

**CAUTION:** To prevent incorrect measurement of the aortic pressure by the guiding catheter, causing an incorrect FFR (Pd/Pa) calculation:

- Flush any contrast remnants from the guide catheter with heparinized saline.
- The insertion tool must be pulled back out of the hemostatic valve.
- The hemostatic valve must be closed during pressure measurement.

**NOTE:** The equalization value (EQ) is reset when Pa or Pd is zeroed, or when the PressureWire is restarted (turned on/off).

## Recording FFR

1. Confirm that the FFR settings are correct for the current procedure. See “[FFR Settings](#)” on page 2-5 for more information.
2. With the PressureWire in position and equalized, advance it across the area to be examined.

The bottom of the graphical area shows the Pd/Pa waveform, along with the Pd/Pa ratio.



Figure 4-7: Pd/Pa waveforms equalizing

3. Induce hyperemia according to standard cathlab procedures.
4. Click the **Record** button (or on the Table Side Controller press the **Proximal Marker** button).

A recording timer is displayed at the bottom of the screen. Click the **Mark** button (or on the Table Side Controller press the **Distal Marker** button) to mark an instant

## Performing an FFR Procedure

### Recording FFR

---

on the recording for later review. A mark appears as a vertical white line on the recording.



Figure 4-8: Recording

5. Record pressure until steady state maximum hyperemic condition is reached, until the hyperemic effect begins to decrease, or until the physician decides to end the recording. When finished, click the **Stop** button (or on the Table Side Controller press the **Proximal Marker** button) to end the recording.

**NOTE:** During the recording, you can mark by clicking on the **Mark** button, or on the Table Side Controller by pressing the **Distal Marker** button

**NOTE:** Turn off the PressureWire transmitter when it is no longer needed and before disposal.

**CAUTION:** When the procedure is complete, handle the PressureWire and all other disposables used during the procedure as potential biohazards. Dispose of in accordance with accepted medical practice and all applicable laws and regulations.


### Reviewing an FFR Recording

When the FFR recording opens, the entire recording is shown. By default, when the recording is opened, the cursor is positioned at the point of the lowest Pd/Pa ratio.

**WARNING:** The system may place the point of FFR at the wrong location due to abnormal heartbeats, for example artifacts in AO (Pa) caused by flushing of guiding catheter or valve opening/closing. The physician should always confirm that the point selected by the system is a valid point for the calculation of FFR.

**CAUTION:** If you have zoomed in on a section of the recording, portions of the waveform are not displayed (the recording extends off screen). The physician should always review the entire unzoomed recording before selecting the point for the calculation of FFR.

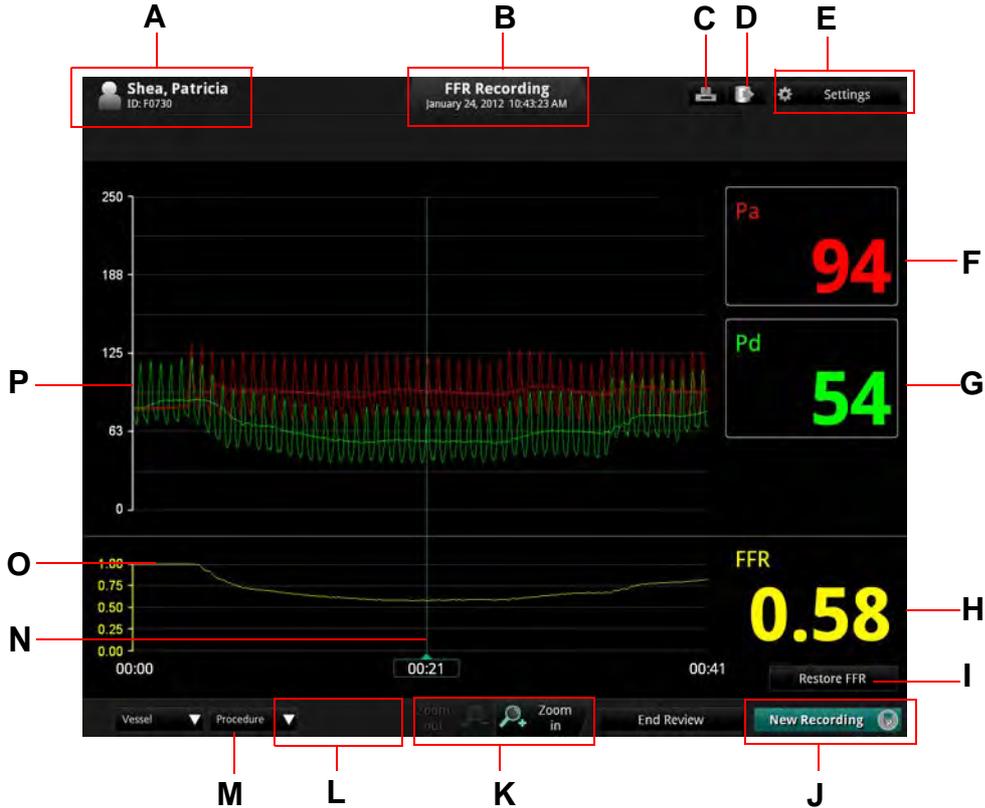Table 4-1: FFR Review Screen



## Performing an FFR Procedure

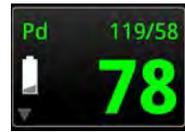
### Reviewing an FFR Recording

---

- A** Patient name and ID.  
See [“Entering New Patient Information” on page 3-6](#) for more information.
- B** Recording date and time.
- C** **Print file to USB** button: Available when a USB drive is connected. Click to print the FFR recording file to a USB drive.
- D** **Export** button: Click to open the **Export Wizard**.  
See [Chapter 8 “Exporting, Importing, and Managing Files”](#) for more information.
- E** **Settings** button: Click to open the FFR Settings menu.  
See [“FFR Settings” on page 2-5](#) for more information.
- F** Pa measurement box. Mean Pa value at the cursor is displayed.
- G** Pd measurement box. Mean Pd value at the cursor is displayed.
- H** The FFR value at the cursor.
- I** **Restore FFR** button: After you have moved the FFR cursor during review; click on **Restore FFR** to reset it to the location it was at when you opened the recording.  
**NOTE:** Upon ending the FFR review, the marker will stay saved to the last point it was on. If it is opened again, the reset button will be grayed out, until the FFR marker is then again moved.
- J** **End Review / New Recording** : Click the **End Review** button to close this window and return to the **Patient Summary** menu.  
**NOTE:** While the system is connected to a PressureWire, the button **Record** appears here. Click the **New Recording** button to close this review and begin a new FFR recording.  
**CAUTION: When closing the recording, the current cursor position and corresponding FFR value is saved. When the recording is re-opened, the cursor appears at the saved position. The cursor may then be moved, and the displayed Pa and Pd pressures and FFR value change to reflect the new cursor position.**
- K** **Zoom in/Zoom out** of the recording, centered on the cursor.
- L** **Procedure** list: Click to open a drop-down list of procedures to describe this recording.
- M** **Vessel** list: Click to open a drop-down list of vessels to describe this recording.
- N** Move the cursor to read the Pa, Pd, and FFR value at any point in the recording.
- O** FFR waveform
- P** Pa and Pd pressure waveforms
-

## PressureWire Troubleshooting

In case of low power in the PressureWire, the light on the PressureWire unit flashes yellow, and a low battery indicator appears in the Pd measurement box.



If there is a problem with your PressureWire, refer to the *PressureWire Instructions for Use* for details.



# Performing an OCT Procedure

# 5

## Overview

The OCT procedure requires two operators; a sterile operator and a non-sterile operator. All steps requiring contact with the Dragonfly Imaging Catheter or the outside of the sterile DOC cover must be performed by the sterile operator. All steps performed within the sterile DOC cover or in direct contact with a keyboard or mouse must be performed by the non-sterile operator.

## Required Material and Equipment

- OPTIS INTEGRATED System
- C7 Dragonfly Imaging Catheter, Dragonfly Duo Imaging Catheter, or Dragonfly OPTIS Imaging Catheter
- Sterile DOC Cover
- 3 mL purge syringe
- Contrast media indicated for coronary use, for purging and flush (allow 15 mL for each run planned)
- 0.014 inch guidewire (with torque device if desired)
- Guide catheter (6 French, 0.068 inch ID or larger, with no side holes)
- Sheath introducer (to match guide catheter)
- Hemostatic Y-Adapter/Connector
- Heparinized, physiologic saline solution, for hydrophilic catheter preparation
- Power injector pump for coronary angiography or manual syringe (capable of injecting 4.0 mL /sec for a total of 14 mL in 3.5 seconds)

## Performing an OCT Procedure Overview

---

### OCT Imaging Overview

1. **Position** - Locate the Dragonfly imaging catheter relative to the target lesion/stent.
2. **Purge** - Clear blood from the catheter lumen, if present, using the attached 3 mL syringe.
3. **Puff** - Inject a small amount (~ 4 mL) of contrast through the guide catheter to evaluate clearance.

If clarity is marginal, check the orientation of the guide catheter and target vessel.

4. **Pullback** - From Live View, select **Enable** to start the imaging process.

**NOTE:** In the left coronary system, guide catheter placement and orientation is key to achieving good contrast flow. This is particularly true in the LCX.

## OCT Operating Modes

During acquisition, the system is divided into Standby View and Live View.

- **Standby View** - The DOC is not rotating the imaging catheter. The last image viewed through the catheter lens is shown on the display.
- **Live View** - The DOC is rotating the imaging catheter at low-speed, and is transmitting images from the catheter lens to the display.

After the image has been recorded, see [Chapter 6 “Reviewing OCT Recordings”](#) for more information on reviewing OCT recordings.

## OCT Recording Types

The system can make two types of recordings: **Pullback (Hi-Res or Survey)** and **Stationary**. See [“Setup - Acquisition Menu” on page 10-2](#) for more information.

**NOTE:** **Hi-Res** is a 54 mm pullback recording with double the frame density as compared with **Survey**. **Survey** is a 75 mm pullback recording with standard frame density.

**NOTE:** The instructions in this manual are for **Pullback** recordings. Differences for **Stationary** recordings are noted where applicable.

## OCT Trigger Types

**Automatic** - The default setting in which the system triggers a pullback automatically when a brief sequence of clear image frames are detected as a result of contrast injection. If the flush injection is not detected within 15 seconds after being enabled, the system returns to **Live View**.

**Manual** - The system does not perform a pullback until the operator clicks the Start Pull-back button on the screen, or presses the **Enable** button on the DOC. If the flush injection is not detected within 15 seconds after being enabled, the system returns to **Live View**.

## Angio-Coregistration

Angio Co-Registration allows the user to align the angiography and OCT recording so that OCT frames correlate with the corresponding angiography position. This functionality allows the user to more easily determine their position in the OCT relative to the lesion and assists in stent location and placement.

## Setting up the OPTIS INTEGRATED System

1. Position the system for use. See [“Positioning the System”](#) on page 2-1.
2. Turn on the system. See [“Power On System Cabinet”](#) on page 2-3.

## Setting up the DOC

The DOC must be readied for use in the sterile environment.

**CAUTION:** Protect the exposed connector inside the DOC from fluids at all times. Fluid contact can disable the DOC and require service.

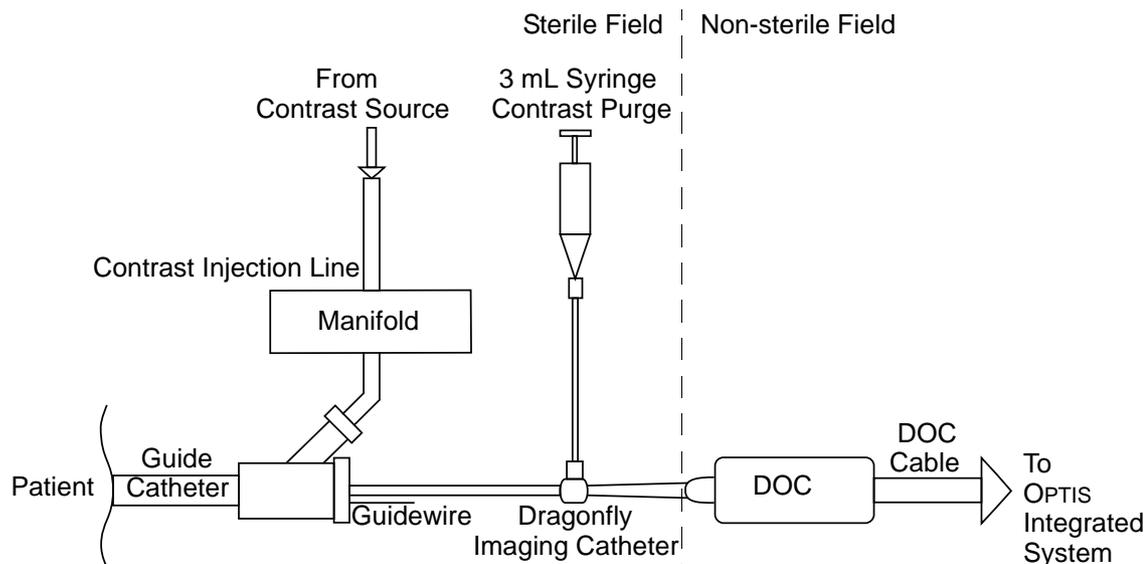


Figure 5-1: DOC Connections

### Non-sterile Operator

1. With assistance from the sterile operator, place the DOC inside a sterile DOC cover.
2. Ensure the cover is extended over the DOC cable to its full length.

### Sterile Operator

3. Connect the flush media (syringe or automated injector) to one port of the guide catheter manifold (see [Figure 5-1](#)).
4. Use the following settings to prepare for flush delivery: (This step can also be performed by Non-Sterile Operator)
  - 4 ml/sec or less flush rate.
  - 14 ml or less total flush volume.
  - If using an automated injector, pressure limit 300 psi, or the nearest available setting.
5. Purge all air from the tubing and manifold following standard practice.

## Setting up the Dragonfly Imaging Catheter

**WARNING:** Prior to use, please review the *Instructions for Use* supplied with the Dragonfly Imaging Catheter for more information.

1. Inspect the packaging of the catheter for damage. Ensure that the catheter has not passed its expiration date.

**CAUTION:** Do not use the catheter if the sterile packaging is compromised.

2. Using sterile technique, carefully remove the Dragonfly catheter from its sterile package.
3. Carefully remove the catheter from the plastic hoop and examine it for visible damage or defects.

**CAUTION:** Do not use the catheter if there are any signs of damage.

4. Moisten the distal segment of the catheter from the tip to approximately 100 cm proximally using heparinized saline to ensure optimal performance of the hydrophilic coating.

**CAUTION:** Use heparinized saline only.

5. Remove the cap from the sidearm luer and attach a 3 ml syringe filled with 100% contrast.

**CAUTION:** Take care in handling the Dragonfly to prevent breaking the fiber-optics within the catheter. Kinking and bending of the catheter can cause damage. While connecting, ensure the proximal catheter segment is straight and aligned with the DOC. Never attempt to connect and operate the catheter while the catheter remains coiled within the hoop.

## **Preparing to Acquire OCT Recordings**

New OCT recordings are started from the **Patient Summary** menu.

1. If necessary, enter the patient's information. See [“Entering New Patient Information” on page 3-6](#).

**NOTE:** Patient information should be entered into the system and selected for use before beginning a recording.

If you are entering a patient through the **Add New Patient** menu, click the **New OCT Recording** button (located at the bottom of **Step 2** of the **Add New Patient** wizard), and continue with [Step 4](#) below.

2. In the **Select Patient** menu, click on the patient.

The **Patient Summary** menu for that patient opens ([Figure 3-3 on page 3-4](#)).

3. In the **Patient Summary** menu, click the **New OCT Recording** button.

**NOTE:** If necessary, enter **Physician name** and **Accession number** for this case, and click **OK**.

4. If **Flush Medium** was set to anything other than **Contrast** for the previous patient, the setting is changed back to **Contrast** and an alert appears.

The system will ask you to confirm or choose another flush medium, then click **OK** to close the alert. If you need to change the Flush Medium setting, see [“Confirm Recording Settings” on page 5-12](#).

5. The screen displays the guidance message “Purge catheter with contrast” (Figure 5-2).



Figure 5-2: Purge Catheter guidance message

**NOTE:** To prevent this guidance message from appearing, check the box next to the message “Do not show this screen again.”

6. Purge the Dragonfly Imaging Catheter's lumen with 100% contrast media from the 3 ml syringe to remove all air from the catheter. Flush until 3-5 drops exit from the catheter's distal tip. Do not remove the syringe from the catheter purge port after flushing.

**CAUTION:** The catheter must be purged prior to connection to the DOC.

## Performing an OCT Procedure Preparing to Acquire OCT Recordings

---

7. Once the 100% contrast media has been injected, click the **Next** button.

The screen displays the guidance message “*Plug catheter into DOC*” (Figure 5-3).



Figure 5-3: Plug Catheter into DOC Guidance Message

8. Insert the hub of the catheter into the port of the DOC, and twist the hub clockwise until secure ( $\frac{1}{8}$  turn).

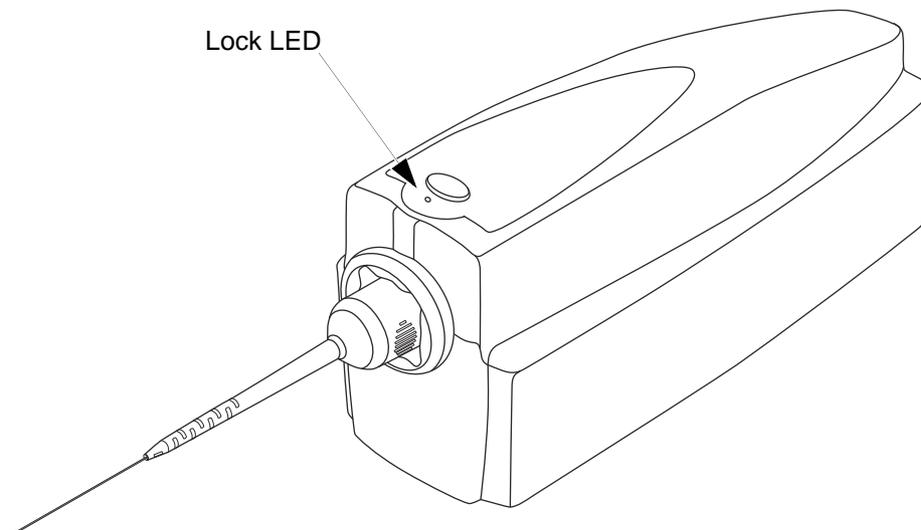


Figure 5-4: Dragonfly Catheter Connected to the DOC

**CAUTION:** Do not insert or remove the catheter while the DOC is scanning. Do not attempt to disconnect the catheter from the DOC while the “lock” LED is blinking as it could damage the catheter or the DOC.

The DOC's “lock” LED flashes and lights as the DOC automatically makes the internal optical fiber connection, and the green **Connecting imaging catheter** progress bar at the bottom of the screen completes.

**NOTE:** If the **Connection Failure** alert appears, unload the catheter (press **unload** on the DOC) and try to load the catheter again. If this alert appears again, replace the catheter.

**NOTE:** If the connection or initial calibration is *stopped*, the **Connection Cancelled** alert is displayed. If this alert appears, unload the catheter (press **unload** on the DOC), click **OK**, and try to load the catheter again.

On screen, the size and position of the outer sheath of the catheter is adjusted with respect to the calibration marks (see [Figure 5-5](#)).

When the connection and initial calibration are complete, the confirmation message, **Catheter connected**, appears. The system goes into Standby View.

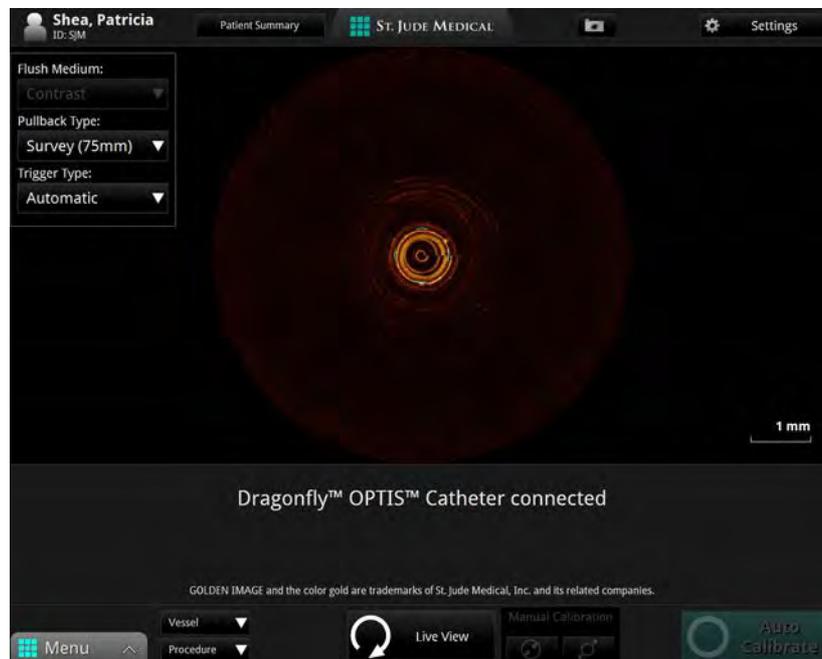


Figure 5-5: Catheter Connected, Initial Calibration done

## Confirm Recording Settings

After the catheter has been connected to the DOC, confirm all settings for this recording.

1. Click the **Settings** button at the top of the screen.

A context-sensitive **Settings** menu opens on the right.

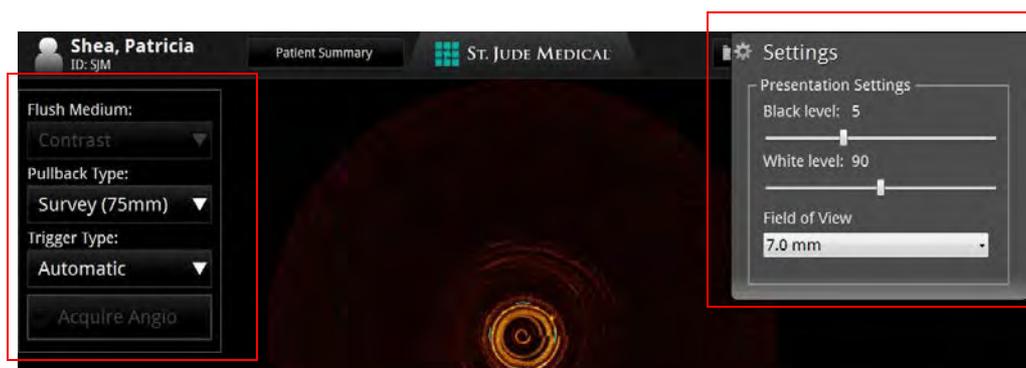


Figure 5-6: OCT Settings Menu (during Recording)

Additional pullback settings are displayed in the menu on the left. For Angio Co-Registration, click on the **Acquire Angio** button. Angiography video input must be present for the **Acquire Angio** button to become live.

2. Confirm all settings are correct for this recording.

**NOTE:** If Flush Medium has been restricted to **Contrast Only**, you cannot change the **Flush Medium** setting. See [“Setup - Acquisition/Other Menu” on page 10-4](#) for more information.

If Recording Type has been set to **Stationary**, you cannot change the **Pullback Length** setting. See [“Setup - Acquisition Menu” on page 10-2](#) for more information.

**CAUTION:** To obtain accurate measurements, be sure the selection for the **Flush Medium** is the same as the medium in which you are imaging.

**NOTE:** Only 100% contrast media is approved for human OCT Imaging.

**NOTE:** Click the **Menu** button at the bottom of the screen to open the OPTIS INTEGRATED System **Setup** dialog box. See [Chapter 10 “User Interface Reference”](#) for more information.

## Dragonfly Imaging Catheter Insertion and Positioning

The guide catheter and guidewire must be inserted into the patient per normal clinical procedures prior to insertion of the Dragonfly Imaging Catheter. The Dragonfly catheter must be properly inserted into the guide catheter to ensure patient safety and proper operation.

Once the Dragonfly catheter has been inserted, its location, and the position of the guide catheter, may be fine-tuned to ensure optimal imaging.

**CAUTION: Ensure that no air is introduced into the system during the imaging catheter insertion.**

1. Ensure that the Dragonfly catheter is not rotating (press **Standby View** if necessary) before the sterile operator begins loading and inserting. If desired, select a **Vessel** and **Procedure** from the fields at the bottom of the screen.
2. Back-load the Dragonfly catheter's rapid-exchange lumen onto the indwelling 0.014 inch guidewire.
3. **POSITION** the catheter according to your catheter's Instructions for Use.

**CAUTION: Observe all advancement and movement of the Dragonfly Imaging Catheter under fluoroscopy. Always advance and withdraw the catheter slowly. Failure to observe device movement fluoroscopically may result in vessel injury or device damage. To ensure proper placement, do not move the guidewire after the Dragonfly catheter is in place.**

**If resistance is encountered during advancement of the Dragonfly Imaging Catheter, stop manipulation and evaluate under fluoroscopy. If the cause of resistance cannot be determined or mitigated, carefully remove the catheter and guidewire as a unit from the patient.**

**Leave the guidewire engaged with the catheter at all times during use. Do not withdraw or advance the guidewire prior to withdrawing the catheter.**

4. Press the **Live View** button to start live-scan imaging (low-speed rotation of the imaging core).

The ring around the **Live View** button on the DOC lights in green.

## Performing an OCT Procedure Dragonfly Imaging Catheter Insertion and Positioning

---

**CAUTION:** Monitor the OCT image for indications of catheter optical failure. If optical failure is suspected, return to Standby View (click the **Standby View** button), remove the catheter, and replace it with a new one.

**NOTE:** A pullback can not be started unless the optical fiber is fully advanced.

5. Calibrate either by pressing the **Auto Calibrate** button or by using the **Manual Calibration** buttons. You can also **Auto Calibrate** by pressing the **Enable** button on the DOC.
6. Verify calibration as follows:

The image is correctly calibrated when the outermost “ring” of the catheter is centered between the 4 calibration marks.

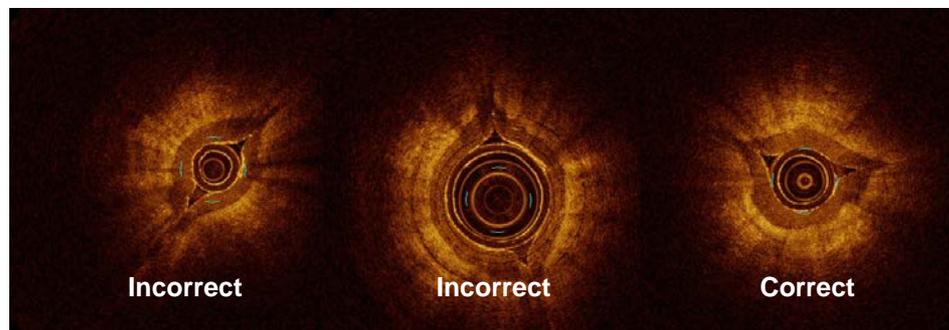


Figure 5-7: Incorrect and Correct Calibration (Dragonfly Duo shown in fingertips)

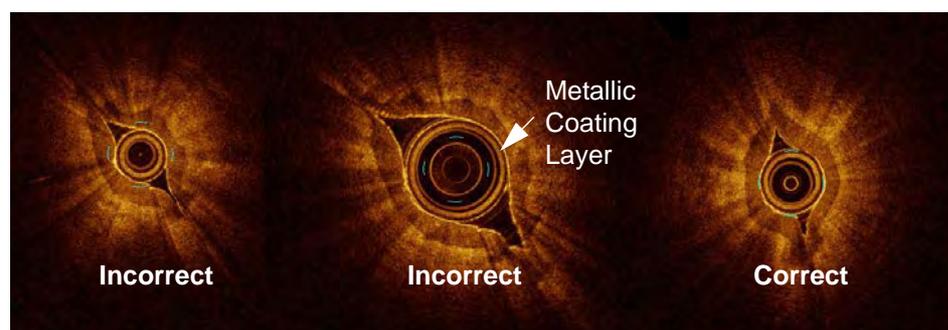


Figure 5-8: Incorrect and Correct Calibration (Dragonfly OPTIS shown in fingertips)

**NOTE:** Ensure accurate calibration prior to a recording. Incorrect calibration may cause early pullback initiation and incorrect measurements.

## Performing an OCT Procedure Dragonfly Imaging Catheter Insertion and Positioning

---

You should adjust the calibration to the best approximate size between the outer diameter of the catheter and the 4 calibration marks.

- To calibrate automatically, press the **Auto-calibrate** button. In the Image Window, the position of the catheter outer sheath is adjusted close to its final correct position with respect to the four calibration marks.
- If necessary, click the Increase/Decrease buttons to calibrate manually. 

7. **PURGE** the Dragonfly Imaging Catheter by injecting ~0.1 ml contrast using the 3 ml syringe to ensure no blood has diffused into the catheter lumen.

**NOTE:** Application of negative pressure to draw blood into the catheter is not recommended. Blood in the catheter lumen will obscure the image and can be difficult to completely purge.

8. Ensure the guide catheter is oriented to preferentially direct contrast flow to the target artery, and verify angiographically that adequate flow of contrast is delivered to the artery.

**NOTE:** **PUFF**-- Injecting a small “puff” of contrast while reviewing the image on screen to verify that adequate flow of contrast is delivered to the artery.

9. The operator can inject the media either by manual flush or automated injector. If using an automated injector, verify the following for flush delivery:
  - 4 ml/sec or less flush rate.
  - 14 ml total flush volume.
  - If using an automated injector, pressure limit 300 psi, or the nearest available setting.
10. Verify the stopcock position on the manifold is set to allow flow from the injection pump into the guide catheter.

## Acquiring Patient Images

The System Display during image acquisition is shown in [Figure 5-9](#).

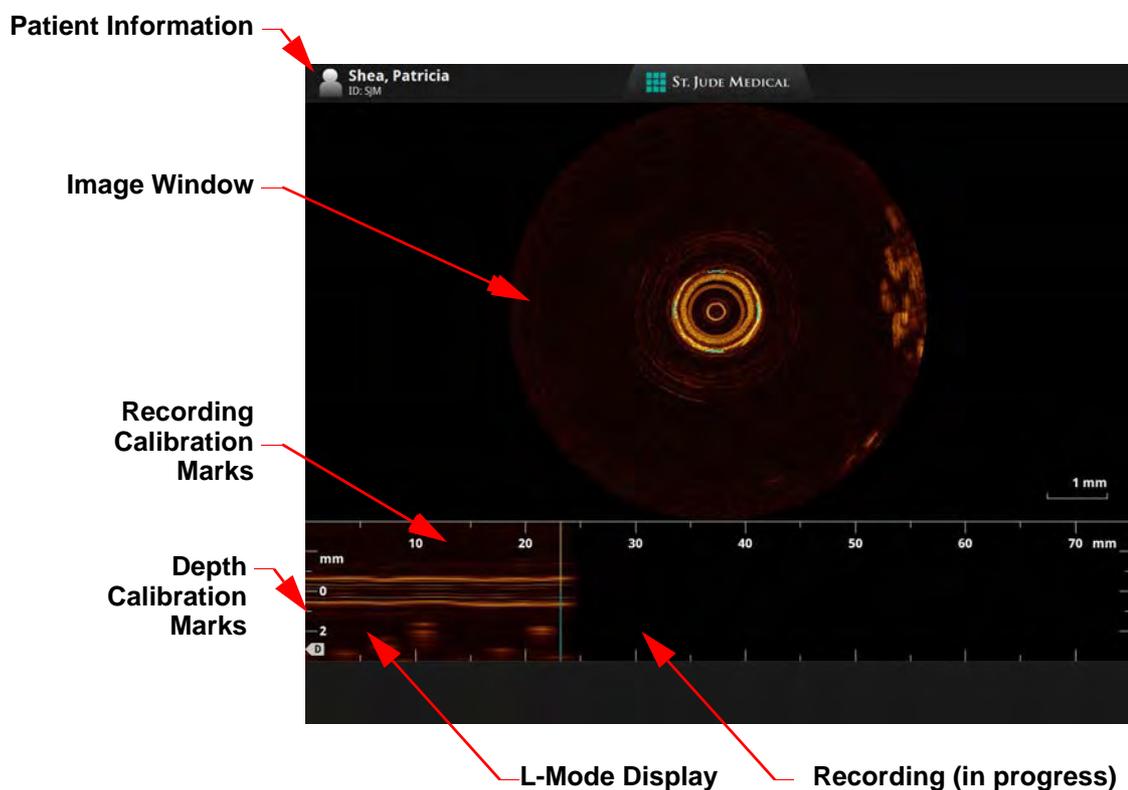


Figure 5-9: System Display - Acquisition

Table 5-1: System Display Description - Acquisition

Patient Information	Displays the Patient name and ID.
Image Window	This display shows the current view by the imaging catheter lens.
Recording Calibration Marks	Recording calibration marks, measured in mm.  Note: In a <b>Stationary</b> recording, this is a timeline, measured in seconds. See <a href="#">“Setup - Acquisition Menu” on page 10-2</a> to set the Recording Type to <b>Stationary</b> .
Depth Calibration Marks	Depth of the scan in mm.
L-Mode display	This displays the approximate lateral appearance of the vessel being scanned.  Note: In a <b>Stationary</b> recording, the “L-Mode” is renamed “Timeline.”
Recording (in progress)	Displays the recording as it is completed.

**NOTE:** If you are performing a **Stationary** recording, the **Enable Pullback** and **Start Pullback** buttons are replaced with a **Start Recording** button. Do not click the **Start Recording** button until you have notified the person injecting the contrast media that the system is ready. Click the **Start Recording** button (or Sterile Operator press the **Enable** button on the DOC) after the contrast is injected. The recording begins immediately.

1. With the system in **Live View**, calibrate either by pressing the **Auto Calibrate** button or by using the **Manual Calibration** buttons. Click the **Enable Pullback** button (or Sterile Operator press the **Enable** button on the DOC) to allow the system to detect initiation of the imaging flush.

You have 15 seconds to initiate the flush for the pullback.

**NOTE:** When acquiring the angiographic images for co-registration with corresponding OCT recording, the imaging equipment must be stationary throughout the cineangiography. For best results, choose the view that reduces occurrences of vessel foreshortening and branch overlapping. This should minimize co-registration inaccuracies.

2. Notify the person injecting the contrast media that the system is enabled.
3. Begin contrast media injection (by manual flush or automated power injection).

The system automatically begins a recording once a brief sequence of clear image frames are detected indicating the area is flushed with contrast media.

**NOTE:** If **Trigger Type** is set to **Manual**, click the **Start Pullback** button or press the **Enable** button on the DOC to initiate pullback.

**NOTE:** If the **Recording Type** is set to **Pullback**, the recording takes approximately 3 seconds. If set to **Stationary**, the recording takes approximately 6 seconds.

**NOTE:** To use Angio Co-Registration with this recording, first click on the **Acquire Angio** button, then step and stay on the **Cine Pedal** while simultaneously capturing the entire pullback recording. Step off the **Cine Pedal** upon completion of the pullback.

If the **Automatically review recordings** option is enabled, the image file is displayed for review.

Review the recording and repeat pullback procedure if needed. If the catheter is still connected to the DOC, the **New Recording** button will be available at the bottom of the screen.

## Removing the Dragonfly Imaging Catheter

All information obtained during imaging is automatically saved with the recording. Once imaging is completed, the data and images may be reviewed (see [Chapter 6 “Reviewing OCT Recordings”](#)) or transferred to removable media for review on a St. Jude Medical Offline Review Workstation (see [Chapter 8 “Exporting, Importing, and Managing Files”](#)).

1. Ensure that the Dragonfly catheter is not rotating (click on the **Standby View** button if necessary) before the sterile operator begins withdrawal and unloading.
2. When all OCT imaging is complete, withdraw the Dragonfly catheter into the guide catheter under fluoroscopic observation.

**CAUTION:** If resistance is encountered during withdrawal of the Dragonfly Imaging Catheter, stop manipulation and evaluate under fluoroscopy. If the cause of resistance cannot be determined or mitigated, carefully remove the catheter and guidewire as a unit from the patient.

3. Remove the Dragonfly catheter from the guide catheter and guidewire.
4. Disconnect the Dragonfly catheter from the DOC by pressing the **Unload** button. When the “lock” LED stops flashing, the catheter is disconnected internally. Twist the catheter hub counter-clockwise to disengage it from the DOC.

**CAUTION:** Handle the Dragonfly Imaging Catheter and all other disposables used during the procedure as potential biohazards. Dispose of in accordance with accepted medical practice and all applicable laws and regulations.

## Troubleshooting OCT Acquisition

### Immediately Stopping DOC Operation

1. Press the **Stop** button on the DOC at any time to immediately stop operation of the DOC. Pressing this button:
  - Turns off the DOC motors, thereby stopping all optical fiber movement, both rotational and longitudinal.
  - Turns off the Imaging Engine light source.
  - Stops recording and freezes the image. If **Automatically review recordings** is checked in the **Configure** tab of the **Setup** dialog box, it automatically plays back the recording once it has been saved.
  - Prevents the imaging catheter lens from automatically advancing to its original distal position.
2. The system computer remains on and the system is available for new operations.

**CAUTION:** If the catheter was stopped before returning to its fully advanced position, the system displays an alert: **Catheter Must Be Advanced to Continue**. Click the **Advance Catheter** button to advance the optical fiber within the catheter. You must advance the catheter before making another OCT recording.

## **Catheter Failure**

**NOTE:** These events can only occur if a catheter is connected to the DOC.

In the event of an imaging catheter failure (complete break of the optical fiber) during either pullback or advance, the system stops automatically (and stops pullback/advance) and displays the following message:

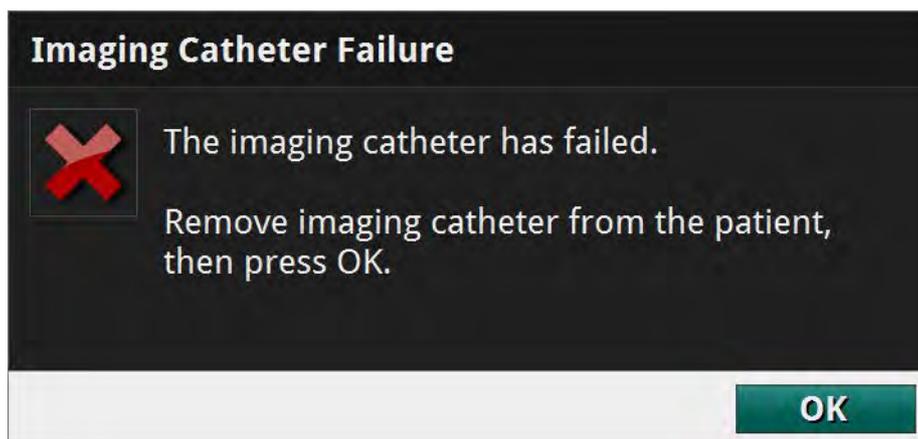


Figure 5-10: Catheter Failure message

**CAUTION:** Do not click the OK button until after you have removed the Dragonfly Imaging Catheter from the patient.

1. Withdraw the Dragonfly into the guide catheter under fluoroscopic observation.

**CAUTION:** If resistance is encountered during withdrawal of the Dragonfly catheter, stop manipulation and evaluate under fluoroscopy. If the cause of resistance cannot be determined or mitigated, carefully remove the catheter and guidewire as a unit from the patient.

2. Remove the Dragonfly catheter from the guide catheter and guidewire.
3. Acknowledge the alert by clicking on **OK**.

After acknowledging the alert, the Safe Unload Guidance displays as shown below.



Figure 5-11: Safe Unload Guidance, Screen 1

4. Turn the catheter hub 1/4-turn counter-clockwise as directed. While the DOC is unloading the catheter, the Guidance displays the screen below:



Figure 5-12: Safe Unload Guidance, Screen 2

## Performing an OCT Procedure Troubleshooting OCT Acquisition

---

5. When the DOC is finished unloading the catheter, the Guidance displays the screen below:

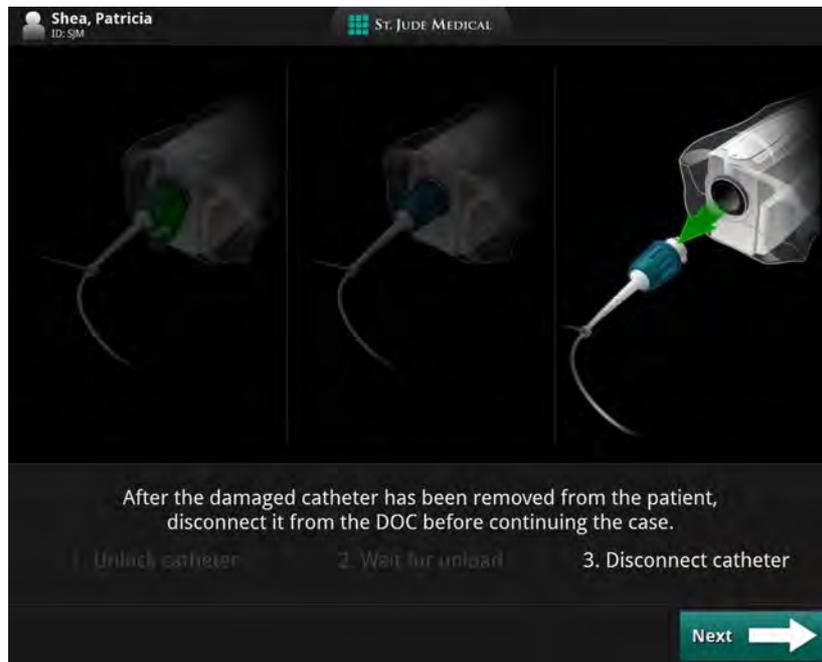


Figure 5-13: Safe Unload Guidance, Screen 3

6. Remove the existing Dragonfly Imaging Catheter in the normal manner. Click **Next**.
7. Load a new Dragonfly Imaging Catheter in the normal manner.

# Reviewing OCT Recordings

---

6

**NOTE:** Making measurements, calculations, and text annotations is covered in [Chapter 7 “Measurements and Annotations”](#). Exporting images is covered in [Chapter 8 “Exporting, Importing, and Managing Files”](#).

## Reviewing OCT Recordings Image Window

### Image Window

When reviewing recordings, the Image Window shows a cross-section of the pullback or still frame ([Table 6-1](#)).

Table 6-1: OCT Display Overview



- A** Patient name and ID.  
See [“Entering New Patient Information”](#) on page 3-6 for more information.
- B** Recording date and time.
- C** **Capture** button: Available on a still frame or paused recording. Click to save the current frame.
- D** **Print** button: Available when a USB drive is connected and the system is displaying a still frame or paused recording. Click to print the current frame to file.  
See [“Printing Still Images”](#) on page 6-13 for more information.
- E** **Export** button: Click to open the **Export Wizard**.  
See [Chapter 8 “Exporting, Importing, and Managing Files”](#) for more information.
- F** **Settings** button: Click to open the (Playback) **Settings** menu.  
See [“Adjust Playback Settings”](#) on page 6-9 for more information.
- G** Frame number: Only visible on a paused recording when Tool Panel is closed.

- H** Tool Panel containing Measurement and Annotation tools: Use these to add measurements, calculations, and add text to recordings and still images.  
See [Chapter 7 “Measurements and Annotations”](#) for more information.
- J** Bookmark controls: Add bookmarks to the L-Mode view.  
See [“Bookmark Controls” on page 6-10](#) for more information.
- K** **End Review / New Recording** : Click the **End Review** button to close this window and return to the **Patient Summary** menu. See [“Patient Summary Menu” on page 3-4](#) for more information.
- I** **NOTE:** While the OPTIS INTEGRATED System is connected to a Dragonfly imaging catheter, the button **New Recording** appears here. Click the **New Recording** button to close this review and begin a new OCT recording.
- L** Playback controls: Control the playback of the OCT recording. Not available with still images.  
See [“Playback Controls” on page 6-6](#) for more information.
- M** **Procedure** list: Click to open a drop-down list of procedures to describe this recording.
- N** **Vessel** list: Click to open a drop-down list of vessels to describe this recording.
- O** **Menu**: Displays the context-sensitive menu. Click to access the **Setup** and playback **Calibration** controls. See [Chapter 10 “User Interface Reference”](#) for more information on the **Setup** Menu; see [“Calibration Adjustment” on page 6-7](#) for more information on the **Calibration** Menu.
- P** L-Mode view: An approximate lateral representation of the vessel for this recording. Not available with still images.  
See [“L-Mode View” on page 6-4](#) for more information.
- Q** Image Window: A cross-section view of the vessel.
- R** View Panel: Control display features including L-Mode and optional features such as 3D and Lumen Profile.
- S** Measurements Panel: This section lists measurements from the current image; click on a measurement to highlight it in the Image Window.
-

## Reviewing OCT Recordings

### L-Mode View


---

## L-Mode View

**NOTE:** If reviewing a **Stationary** recording, the L-Mode display is renamed *Timeline*. The Timeline represents the stationary view of the catheter during the six seconds of recording.

During the OCT pullback recording, the system captures evenly spaced cross-section images and uses them to construct a lateral view of the vessel anatomy. The lateral view is shown in the L-Mode display in the lower portion of the screen; the distal portion of the recording is to the left.

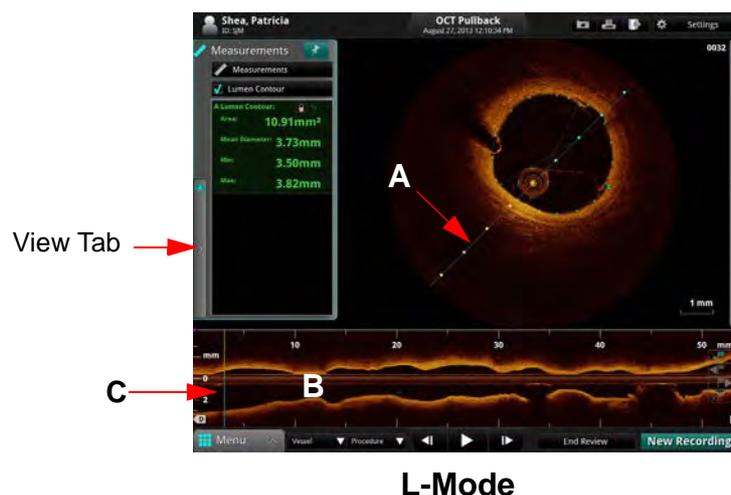Table 6-2: L-Mode view



- A** Cut-Plane Indicator.
- B** L-Mode display.
- C** Current Frame Indicator. Click and drag this to change the lateral view shown in the L-Mode display.

**NOTE:** The width of the L-Mode window represents the entire range of the recording.

## Limitations of L-Mode Data


Please be aware of the following limitations of L-Mode data:

- Due to the catheter's arbitrary position within the vessel, L-Mode data may not be representative of the actual vessel lumen.
- The vessel diameter shown in the L-Mode reconstruction may appear significantly smaller than the actual diameter when the catheter position is off center and close to the vessel wall. To avoid misinterpreting the image when this occurs, rotate the cut-plane to examine all views of the L-Mode.
- Although a vessel may curve, the L-Mode view always appears straight because of the limitations of reconstruction.
- Artifacts caused by the relative motion of the catheter and the vessel often result in a saw-toothed appearance of the reconstruction and can lead to misinterpretation by inexperienced users.
- Shortening or lengthening artifacts in the L-Mode reconstruction may occur due to the relative motion of the Dragonfly catheter with respect to the coronary artery caused by the patient's heart motion.

## Playback Controls

Table 6-3: Playback Controls

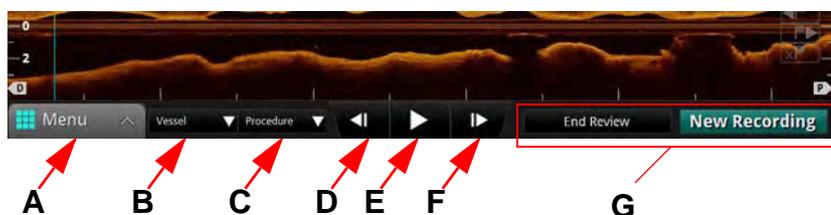---



- A** **Menu**: Displays the context-sensitive menu. Click to access the Setup and playback Calibration controls. Only available with a still image or paused recording.
- B** **Vessel list**: Click to open a drop-down list of vessels to describe this recording.
- C** **Procedure list**: Click to open a drop-down list of procedures to describe this recording.
- D** **Step Backward** : Click to move the recording back one frame at a time. Click and hold this button to move back rapidly. Not available with still images.
- E** **Play/Pause** : Plays or Pauses the recording. If you pause a recording, the frame number is displayed in the upper right corner of the image area. Not available with still images.
- F** **Step Forward** : Click to move the recording forward one frame at a time. Click and hold this button to move forward rapidly. Not available with still images.
- G** **End Review / New Recording** : Click the **End Review** button to close this window and return to the **Patient Summary** menu.

**NOTE:** While the system is connected to a Dragonfly Imaging Catheter, the button **New Recording** appears here. Click the **New Recording** button to close this review and begin a new OCT recording.

---

## Calibration Adjustment

Use the **Adjust Calibration** tool to adjust the calibration of the recording.

1. Use the Playback Controls at the bottom of the screen to pause the playback.
2. Hover or click on the **Menu** button at the bottom of the screen and select the **Calibration** option.

The **Calibration** tool opens. The image is zoomed and the calibration circle with two control points is sized approximate to the outside diameter of the catheter.

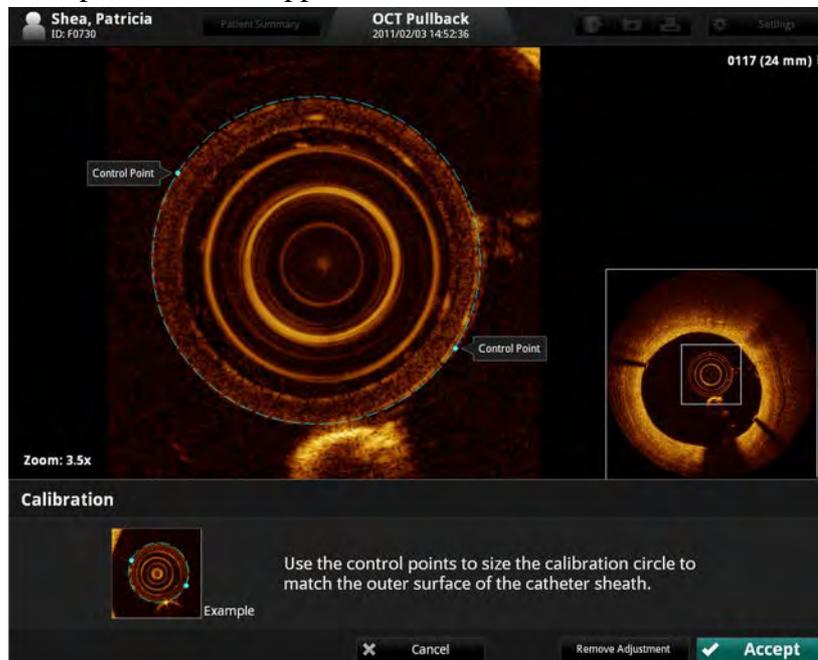


Figure 6-1: Playback Calibration (DragonFly Duo shown, in progress)

3. Click and drag on the control points until the calibration circle traces the outside diameter of the catheter as shown in the sample picture at the bottom of the screen.

**NOTE:** Adjustments made here are for *size*, not *alignment*. The catheter may not appear centered in all frames during playback. This is normal.

4. Click the **Accept** button save the calibration, or click **Cancel** to close the calibration tool without saving the adjustment.

## Reviewing OCT Recordings Calibration Adjustment

---

**NOTE:** If necessary, you can return the recording to its default calibration by clicking the **Remove Adjustment** button. If you remove or change the calibration, any changes that you have made to the automatically-generated Lumen Profile contours are reset. See [“Lumen Profile Display Option” on page 7-20](#) for more information.

## Adjust Playback Settings

1. Click the **Settings** button at the top of the screen.  
  
A context-sensitive menu opens.
2. Click and drag the **Playback Speed** slider bar to set the playback speed.
3. Click and drag the **Rotation** slider bar to rotate the image shown in the cross-section view.  
  
**NOTE:** When an image is rotated, the L-Mode cut-plane marker (if visible) and any displayed measurements and annotations are also rotated.
4. Click and drag the **Black level** and **White level** slider bars to set the black and white balance in the image.
5. Click the arrow on the **Field of View** drop-down menu to display the list of diameter sizes. The default setting is **7.0 mm**.

The **Field of View** setting controls the size of the image displayed on screen. A smaller Field of View setting equals a larger magnification. Click a size to select it.

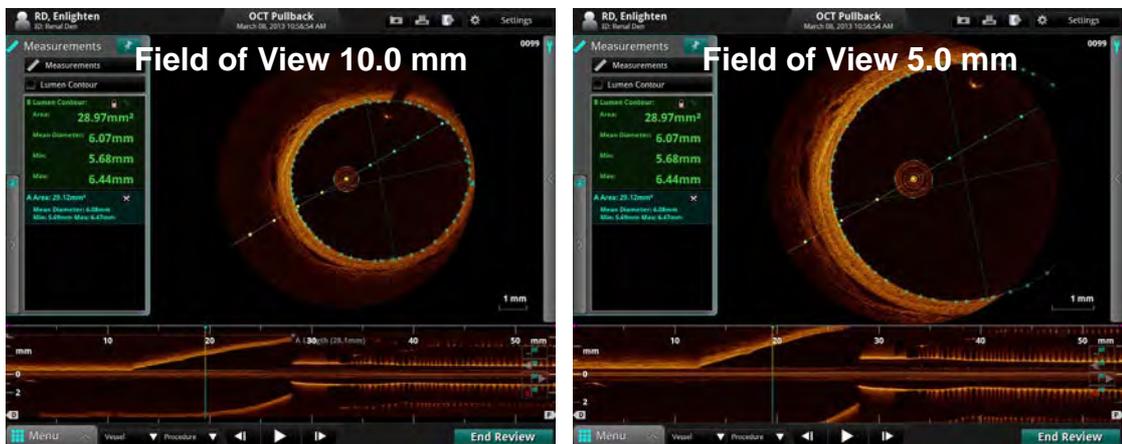


Figure 6-2: Field of View Settings

**NOTE:** Clicking on the image will also change the zoom level. See [“Quick Zoom” on page 7-16](#).


## Bookmark Controls

The System allows you to add bookmarks to mark frames for further review.


**NOTE:** Bookmark controls are only available when an L-Mode view is displayed in a recording.

Table 6-4: Bookmark Controls

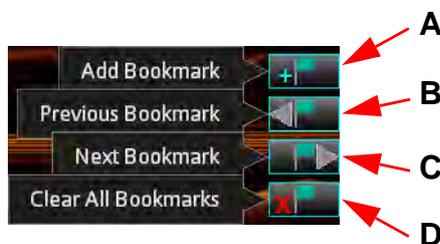---



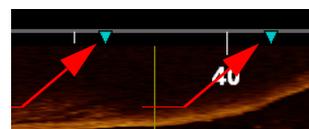**A Add/Remove Bookmark :** Applies or removes a bookmark from the current frame. After a bookmark has been applied to a frame, the position of the bookmark is indicated in the L-Mode view with a green triangle.



Click on a bookmark indicator to jump to the bookmarked frame.

**NOTE:** Frames with measurements and annotations are bookmarked automatically. See [Chapter 7 “Measurements and Annotations”](#) for more information.

- B Previous Bookmark :** Seeks backward to the previous bookmarked frame. If there are no previous bookmarks, it continues seeking from the end of the recording. Unavailable if there are no bookmarks.
- C Next Bookmark :** Seeks forward to the next bookmarked frame. If there are no subsequent bookmarks, it continues seeking from the beginning of the recording. Unavailable if there are no bookmarks.
- D Clear All Bookmarks :** Clears all bookmarks from the current recording. Unavailable if there are no bookmarks.
-

## Setting Playback Range

By default, the system plays the entire length of a pullback during the review. You can shorten the length of the playback by moving the ends of the playback range.

**NOTE:** Playback range markers are only available on recordings in the L-Mode view. By default, they are at the distal and proximal ends of the recording.

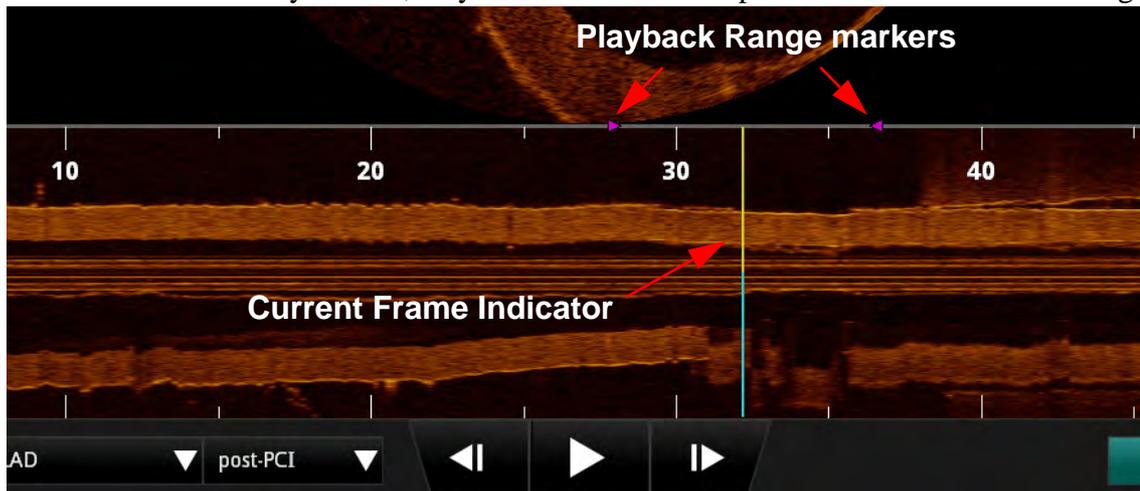


Figure 6-3: Adjusted Playback Range

To change the playback range:

1. Use the Playback Controls at the bottom of the screen to pause the playback.
2. Click on the range markers (purple triangles) and position them before and after the area to be played during review.

**NOTE:** If the recording is exported as an AVI file, this range determines the length of the exported recording.

**NOTE:** If the L-Mode view is changed (for example, turned **Off**), the playback range is restored to the full length of the recording.

## Exporting a Recording or Still Frame



Refer to [Chapter 8 “Exporting, Importing, and Managing Files”](#) for more information on exporting a recording or still frame.

## Capturing Still Images



You can use the **Capture** button to save a still image from a recording, or to save a copy of an existing still frame. All measurements and annotations on that screen are saved with the captured image.

**NOTE:** The **Capture** button is unavailable while a recording is playing. Pause the playback at the frame that you want to capture.

## Saving a Still Image

1. Display the frame that you want to capture.
2. Click the **Capture** button.

A confirmation message appears, and the captured image appears in the **Patient Summary** menu.

**NOTE:** The still image is saved with the same date as the source file. It is grouped in the **Patient Summary** menu with other recordings and images of the same date.

In the **Patient Summary** menu, the system adds the frame number to the title of the captured image (for example, “*OCT Frame 145*”).

## Printing Still Images



You can use the **Print** button to print the current screen (including still frame, L-Mode, and any visible annotations) to an attached USB drive. All measurements and annotations on that screen are saved with the captured image.

**NOTE:** The **Print** button is unavailable while a recording is playing. Pause the playback at the frame that you want to print.

**NOTE:** If you are printing the file to a USB drive, you can set the file format. See [“Setup - Print Menu” on page 10-25](#) for more information.

## Printing a Still Image

1. Display the frame that you want to print.
2. Click the **Print** button.

A confirmation message appears, and the file is saved on the attached USB drive.

## Reviewing with Angio-Coregistration

Only pullback images that contain angio-coregistration data can be used for this type of review. Identify the presence of angio-coregistration data by locating a black and white angiographic thumbnail image in the upper-left corner of the pullback on the **Patient Summary** screen.

1. From the **Patient Summary** screen, identify a pullback that contains a black and white angiographic thumbnail image in the upper-left corner.

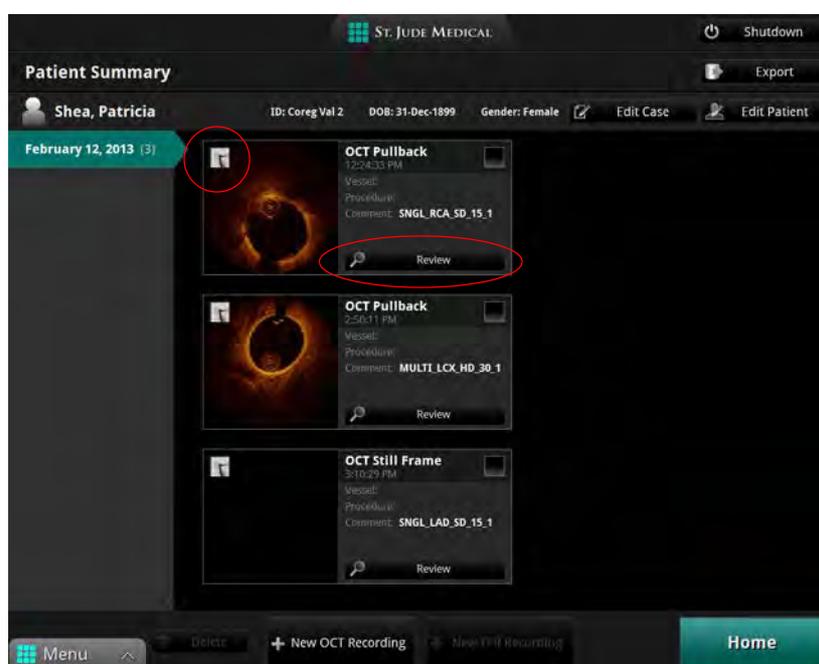


Figure 6-4: Angio-Coregistration Thumbnail Image

2. Click the review button. The OCT review begins to playback.



Figure 6-5: View Menu, Angio Co-Registration Button

3. Pause the recording and click the Angio Co-Registration button under the **View** menu to initiate the Angio Co-Registration process. The angiography displays in the left window.

**NOTE:** To view the angiography along with the MLA data, click the MLA box in the **View** menu. **Mean Diameter** or **Area** is automatically selected when MLA is clicked.

**Reviewing OCT Recordings**  
**Reviewing with Angio-Coregistration**

---

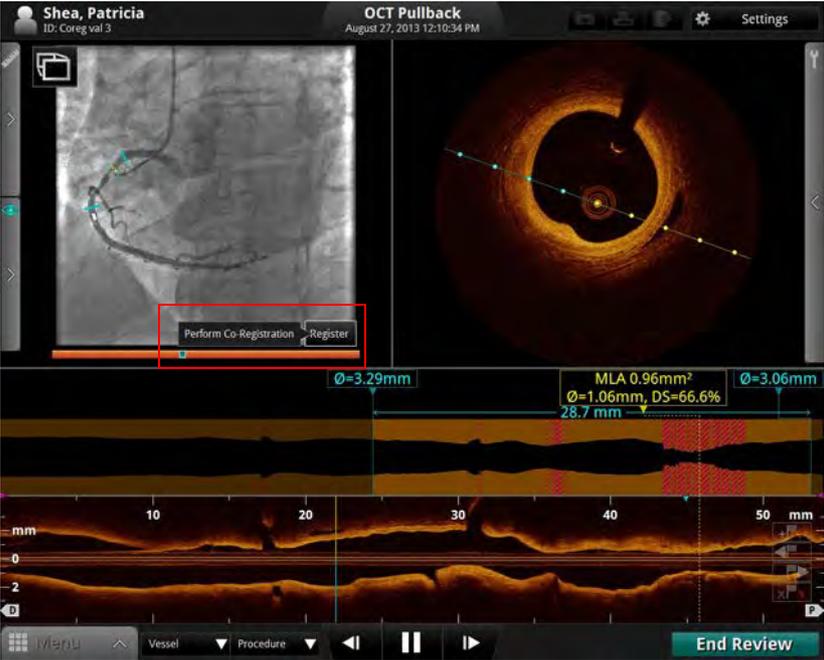


Figure 6-6: Angio Co-Registration: Register

- 4. Click on the **Register** button to initiate the co-registration. Step 1 of the **Angio Co-Registration Guidance** displays.



Figure 6-7: Angio Co-Registration Guidance, Step 1

5. Set the first (distal) control point by placing the cursor in the vessel of interest, on the guidewire (or just proximal of the guidewire), and clicking. A white control point and Step 2 of the Angio Co-Registration Guidance displays.



Figure 6-8: Angio Co-Registration Guidance, Step 2

6. Moving distal to proximal, set at least one additional control point by placing the cursor in the vessel of interest, distally near the guide catheter tip and clicking. Additional white control point(s) display. The system draws a line which traces the vessel of interest and connects the control points. The **Continue** button displays.

## Reviewing OCT Recordings

### Reviewing with Angio-Coregistration

---



Figure 6-9: Angio Co-Registration Guidance, Step 2 with Trace

**NOTE:** If after placing the proximal point, an incorrect vessel path is traced, click the **Undo** button, then replace the proximal control point. Follow the on-screen guidance.

7. Verify that the trace is in the vessel of interest. Click the **Continue** button. Angio Co-Registration Guidance, Step 3 displays.



Figure 6-10: Angio Co-Registration Guidance, Step 3

8. Click the **Confirm** button to accept the path within the vessel of interest. The **Please wait for co-registration** screen displays while the system correlates the angiography with the OCT frames. When the co-registration is complete, the **Coregistration completed successfully** screen displays.

**NOTE:** If either the selected proximal or distal points do not allow the system to generate a full Angio Co-Registration, a **Warning** screen displays. Click **Cancel**. Following the on-screen guidance, replace both distal and proximal control points along the vessel of interest, then click the **Confirm** button.

# Reviewing OCT Recordings

## Reviewing with Angio-Coregistration

---



Figure 6-11: Angio Co-registration completed successfully screen

- 9. Click the **Accept** button to use the successful Angio Co-Registration. The main screen displays the co-registered images.

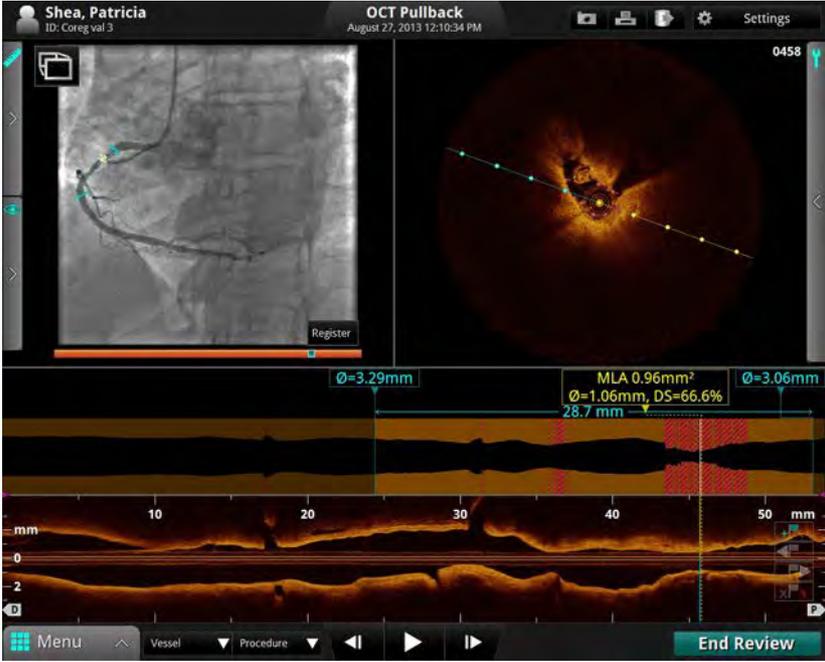


Figure 6-12: Main Screen showing Angio Co-Registration

**NOTE:** To use the Quick Zoom feature (see [“Quick Zoom” on page 7-16](#)) on the Angiography window, simply click on the angiography image.

**Reviewing OCT Recordings**  
**Reviewing with Angio-Coregistration**

---

## Measurements and Text Callouts in the Image Files

Measurements and text callouts that are added to images do not change the underlying image data. The unannotated image is always preserved and can be reviewed without the measurements and callouts. All measurements and annotations added during the procedure are preserved in the file.

Measurements and text callouts can be modified or deleted, and pen color, line width and point size can be set. In addition, pen color can be set to automatically cycle so that subsequent measurements are displayed in different colors. For information on setting pen color, line width, and point size, see [“Setup - Measurements/Labels Menu”](#) on page 10-24.

**CAUTION:** If you want to make measurements on files which will be exported to standard formats, you must make the measurements **BEFORE** exporting the images. Using non-OCT software to measure standard format images will not produce accurate measurements.

**CAUTION:** Do not use images that have been exported to JPEG or Compressed AVI formats for clinical decision making. These formats use compression methods that may degrade the image quality.

**NOTE:** All measurements and calculations can be made in the cross-section view of the OCT image, but only horizontal length measurements and text annotations can be made in the L-Mode view.

## Measurement and Annotation Tools

When the OCT playback is paused, or displaying a still image, the measurement and annotation tools can be accessed on screen by hovering over the **Tools** Panel (wrench symbol). As measurements are added to the image, they are labeled with successive control letters (A, B, C, etc., up to 26 measurements per frame).

**NOTE:** Tools appear black when they are unavailable.

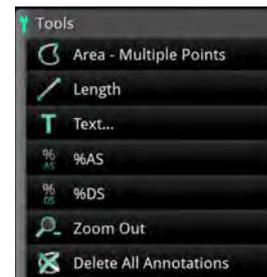


Figure 7-1: Tools

Table 7-1: Measurement and Annotation Tool Functions

---

 **Area - Multiple Points** : Manually place points to trace and measure an area on a frame. This also generates the min., max., and mean diameter measurements for the given area.

 **Length** : Make a distance measurement on a frame.

 **Text...** : Opens the **Enter Note** menu to enter text at the cursor position.

 **%AS (Percent Area Stenosis)** : Calculates the relative size of two areas drawn in the current frame.

**NOTE:** If there are fewer than two areas on the current frame, this button is not available.

 **%DS (Percent Diameter Stenosis)** : Calculates the relative size of two lengths (or diameters) drawn in the current frame.

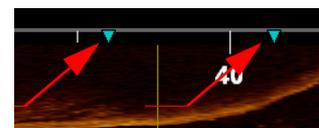
**NOTE:** If there are fewer than two lengths on the current frame, this button is not available.

 **Zoom - / Zoom +** : Switches between zoom levels. See [“Adjust Playback Settings” on page 6-9](#) for more information.

 **Delete All Annotations** : Click to delete all annotations and measurements.

---

**NOTE:** After a measurement or text annotation has been applied to a frame, the position of the bookmark is indicated in the L-Mode view with a green triangle. Click on a bookmark to jump to it.



## Verifying Calibration

Before making measurements, you must verify the calibration. If the calibration for this recording or still image has not been adjusted, or if it has been reset, the calibration tool appears when you begin to take a measurement. See [“Calibration Adjustment” on page 6-7](#) for more information.

Once the calibration is accepted, you can begin taking measurements.

**NOTE:** Calibration before making measurements should not be necessary when using the Dragonfly OPTIS catheter. The OPTIS Integrated software in conjunction with the Dragonfly OPTIS catheter uses *Continuous Calibration* technology to eliminate the need for this step in the workflow.

## Techniques to Improve Measurement Accuracy

Follow these guidelines to improve measurement accuracy:

- Before making measurements, use the Zoom function to zoom into the region of interest until you can clearly see borders and other features you want to measure (see [“Field of View” on page 7-14](#)).
- Place the measurement cursor correctly on the image, using the same measurement technique each time you perform the same type of measurement.
- Avoid making measurements in areas that have artifacts that disguise tissue.

## Measurements and Annotations in the L-Mode View

All annotations and measurements are typically made in the cross-section view area of the Image Window. Only horizontal measurements, useful for determining pullback distance, are permitted in the L-Mode view.

**CAUTION:** Artifacts may result in misrepresentation of L-Mode data, so L-Mode is not recommended for quantization of clinical information.

## Length Measurements

-  The system calculates and displays length as the distance in millimeters (mm) between 2 points placed on an image in either the cross-section or L-Mode views.

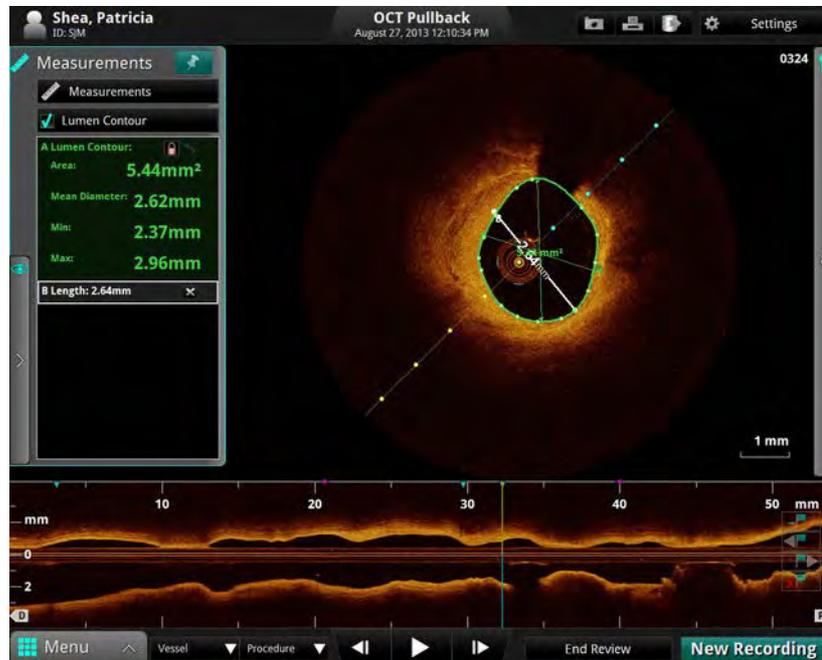


Figure 7-2: Length Measurement

### Making a Length (Distance) Measurement

1. Select the still image or paused recording that you want to measure.
2. Click on the **Length** button (Figure 7-2). The cursor changes to a pen.
3. Click anywhere in the Image Window (cross-section view or L-Mode view) to place the starting point for the length measurement.

**NOTE:** Both start and end points must be in the same view. For example, if the starting point is placed in the cross-section view, the end point must also be in the cross-section view. To cancel the measurement, press the <Esc> key or click the **Cancel** button.



4. Use the mouse to position the pen cursor at the end point and click to set the point.

The completed, labeled distance measurement result appears in the Annotations Panel. If you position the cursor over the line, the length is shown over the center of the line.

## Area Measurements



You can manually create a closed area trace of the lumen contour in the cross-section view. The area is calculated using Green's Theorem, and is displayed in the frame in  $\text{mm}^2$  with the minimum and maximum diameter chords.

**NOTE:** If the Automatic MLA and %DS option is enabled and you are reviewing a pullback recording, your system automatically adds a “Lumen Profile” measurement on each frame.

**NOTE:** The minimum diameter chord has arrowheads that point inward (toward the diameter chord). The maximum diameter chord has arrowheads that point outward.

### Making a Manual Area Measurement

1. Select the still image or paused recording that you want to measure.
2. Click on the **Area - Multiple Points** button (Figure 7-1). The cursor changes to a pen.
3. Use the mouse to position the pen cursor at the desired starting point for the area in the cross-section view and click to set the point.

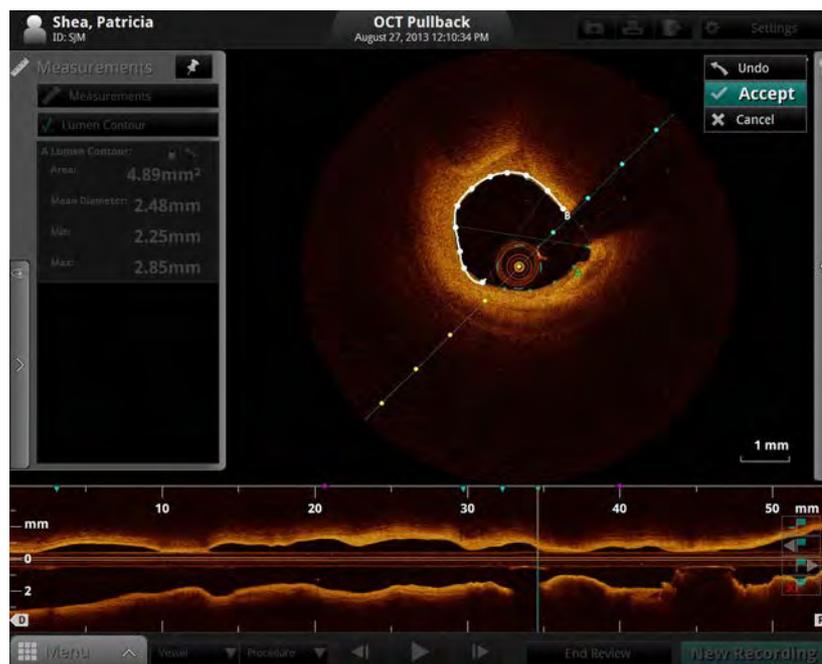


Figure 7-3: Manual Area Measurement (in progress)

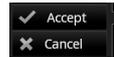
**NOTE:** Click the **Undo** button to the left of the toolbar to delete the most recent point added to the image.



4. Continue to add points with the pen cursor until you have accurately traced the area to be measured. You can place as many points as you like around the border of the desired area. At least three points are required.

**NOTE:** The first two points are connected by a straight line. When you place subsequent points, the straight line becomes a smooth curve connecting all points but not closing the area until you click the **Accept** button, or the last point is sufficiently close to the initial point for the system to automatically complete the area. Placing more points increases measurement accuracy.

**NOTE:** You must click either **Accept** or **Cancel** to complete the measurement. Click **Cancel** to completely erase the measurement.



5. Click the **Accept** button to the left of the toolbar to save the measurement, or click the **Cancel** button to cancel it.

## Adding Text Callouts



You can add text callouts to a single frame or to the entire recording.



Figure 7-4: Text Callouts

## Adding Text Callouts

1. Select the still image or paused recording that you want to annotate with text.
2. Click on the **Text...** button (Figure 7-1).

The cursor changes to a Text marker (“A”).

3. Place the cursor where you want the text to be displayed.

The **Enter Note** menu (see Figure 7-5) appears.



Figure 7-5: Enter Note Dialog Box

4. Type the desired text into the box.
5. If you want to display the text on all frames, click the **Apply to all frames** checkbox. (Not applicable for text added in the L-Mode view.)
6. Click **OK** to approve the note or click **Cancel** to cancel the note.

The text appears, including a callout line beneath the text.

7. To change the position of the text on screen, click and drag the text to the desired position.
8. To change the position of the callout line, click the end of the line and drag it to the desired position.

## The %AS Calculation



The %AS (Percent Area Stenosis) calculation calculates the percentage size of a smaller area with relationship to a larger area, typically the area of the inner border of a vessel compared to the area of the outer border.

### Formula for %AS Calculation

$$\%AS = (\text{First Area} - \text{Second Area}) / \text{First Area} * 100$$

To make this calculation, you must have two areas drawn on the current frame.

### Make a %AS Calculation

1. Make an area measurement.
2. If necessary, make a second area measurement in the frame.
3. Click on the %AS button (Figure 7-1) to display the **Select Area Measurement** dialog box.

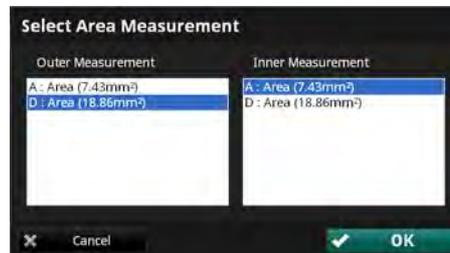


Figure 7-6: Select Area Measurement Dialog Box

4. In the **Outer Measurement** list, click the larger area measurement.
5. In the **Inner Measurement** list, click the smaller area measurement.
6. Click **OK**.

The % Area Stenosis calculation is displayed immediately below the previous calculations in the Annotation Panel as shown in [Figure 7-7](#).

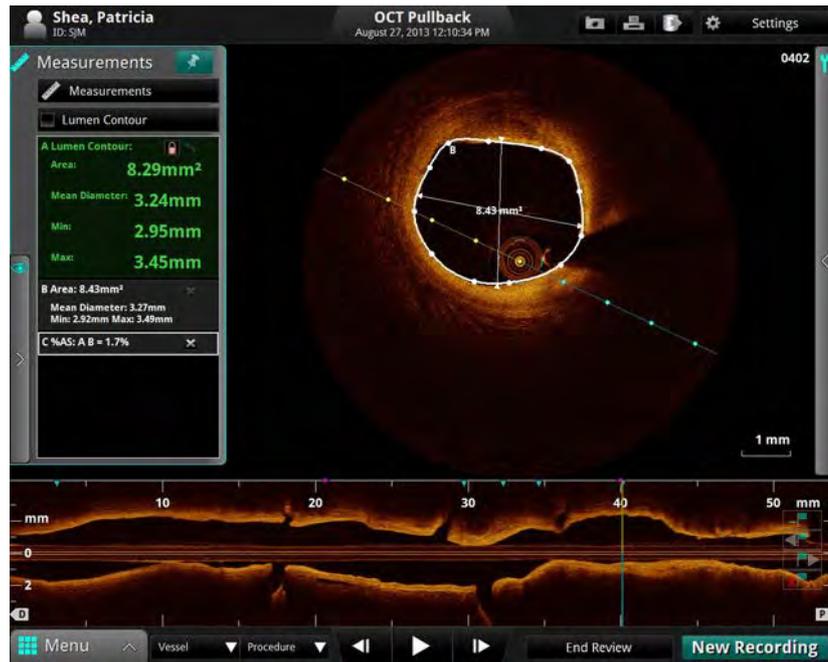


Figure 7-7: %AS Calculation

**NOTE:** If the second measurement you select is larger than the first measurement, a %AS value will not be displayed. Instead, the Annotation Panel displays an error message, similar to the one shown in [Figure 7-8](#).

```
A Area: 8.94mm2
Mean Diameter: 3.37mm
Min: 3.22mm, Max: 3.55mm
B Area: 12.59mm2
Mean Diameter: 3.99mm
Min: 3.72mm, Max: 4.23mm
C %AS N/A: Outer Area (A) < Inner Area (B)
```

Figure 7-8: %AS Error Message

## The %DS Calculation



The %DS (Percent Diameter Stenosis) calculation calculates the percentage size of a smaller length with relationship to a larger one, typically the diameter of the inner border of a vessel compared to the diameter of the outer border.

### Formula for %DS Calculation

$$\%DS = (\text{First Length} - \text{Second Length}) / \text{First Length} * 100$$

To make this calculation you must first draw and specify two lengths, one for the smaller diameter and one for the larger diameter.

### Make a %DS Calculation

1. Make two length measurements in the current frame.
2. Click on the %DS button (Figure 7-1) to display the **Select Length Measurement** dialog box.

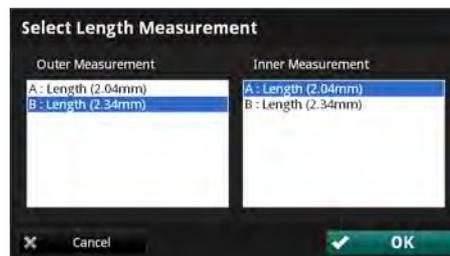


Figure 7-9: Select Length Measurement Dialog Box

3. In the **Outer Measurement** list, click the longer length.
4. In the **Inner Measurement** list, click the shorter length.
5. Click **OK**.

The % Diameter Stenosis calculation is displayed immediately below the previous calculations in the Annotations Panel, as shown in [Figure 7-10](#).

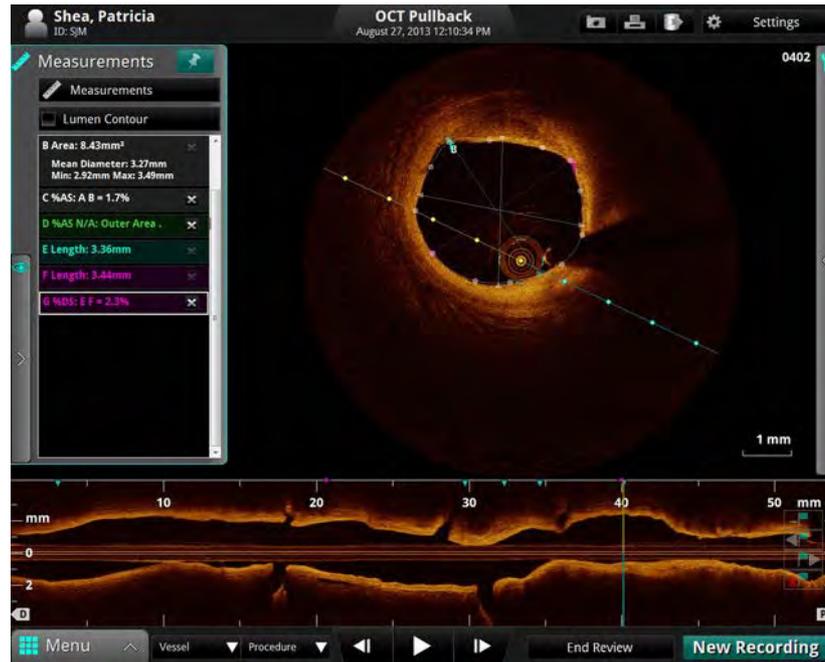


Figure 7-10: %DS Calculation

**NOTE:** If the second measurement is larger than the first measurement you select, a %DS value will not be displayed. Instead, the Annotations Panel displays an error message, similar to the one shown in [Figure 7-8](#) for %AS.

## **Field of View**

You can zoom out on a recording or still image in order to see the full lumen on a large vessel, centered on the catheter.

**NOTE:** Using either Zoom function does not change the image data that is saved; it merely changes the portion of the data that is displayed.

### **Increase/Decrease Field of View**



Click to increase the field of view on the recording or still image. Click again to return to normal image size.

**NOTE:** Depending on the current Flush Medium, the image may not occupy the entire 10.5 mm diameter.

You can change the magnification of the zoomed field. See [“Adjust Playback Settings” on page 6-9](#) for more information.

## Zooming In Manually

You can enlarge any section of the image manually.

**NOTE:** This magnification function cannot be used in the L-Mode display.

1. Place the cursor over the image area you want to enlarge.
2. Click and drag the mouse across an area to magnify it.

When you release the mouse, the system zooms in to the selected rectangle, displaying the magnification factor used (ex. **Zoom: 2.5x**, see [Figure 7-11](#)).

The system adds an overall view of the image to the right.



Figure 7-11: Zooming an Image

3. To return to normal imaging, click anywhere in the zoomed image.

## Quick Zoom



The **Quick Zoom** feature provides two pre-set levels of magnification by simply clicking on the image with the left mouse button. Hover over the cross-sectional view so your cursor changes into a magnifying glass icon. Click once to zoom in to the first level; click a second time to zoom in to the second level; and then click again to return to normal magnification (1.0X).

The pre-set **Quick Zoom Factors** can be customized in the **Settings** menu **Display** option. By default the **First Click** is set to 3.0X. The **Second Click** is set to 6.0X.

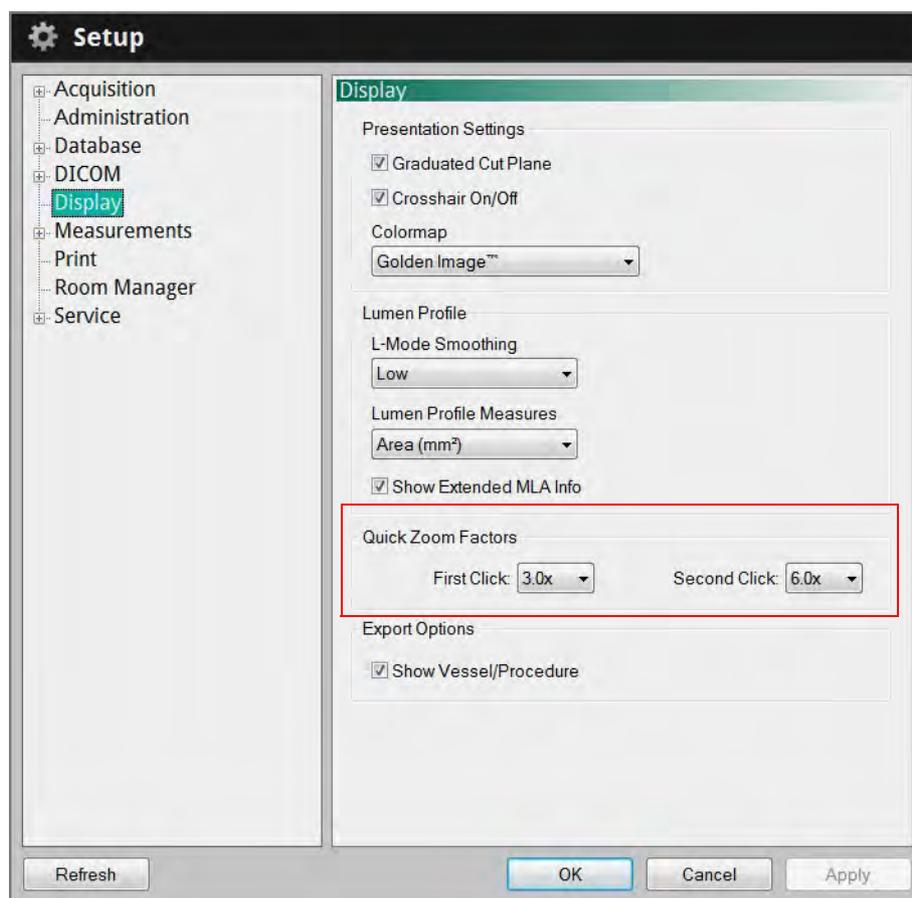


Figure 7-12: Setup Menu, Display Option

## Editing Measurements and Annotations

All measurements and annotations can be moved, deleted, or edited.

- Area and length measurements can be adjusted by moving or deleting control points.
- Text annotations can be moved or deleted.

**NOTE:** When a measurement is changed, calculations that depend on that measurement are updated automatically.

### **Moving Individual Points**

To move an individual point, click on the point and drag it to a new location.

**NOTE:** When you select a point, the point changes from a circle to a square, indicating that it may now be moved.

**NOTE:** If you move a point on an automatically-generated Lumen Profile contour, an arrow icon appears in the corner of the measurement listed in the Annotations Panel. Click the arrow to reset the contour on that frame back to the automatically-generated Lumen Profile contour. See [“Lumen Profile Display Option” on page 7-20](#) for more information. 

### **Adding Points to a Multiple Point Area**

To add additional points to a multiple point area, click on the curve between two points. A new point will be added where you click.

**NOTE:** When you position the mouse over a location where a point can be added, the cursor changes from an arrow to a hand.

### **Deleting Points from a Multiple Point Area**

To delete selected points in a multiple point area, click on the desired point and press the <Del> key on the keyboard.

**NOTE:** When you select a point, it changes from a circle to a square, indicating that it may now be deleted.

## Deleting Individual Measurements or Text Callouts

**NOTE:** If a measurement is in use with a calculation (for example, an area used by **%AS**), the individual measurement cannot be deleted until the calculation is deleted.

To delete a measurement:

- Click on the “x” next to it in the Annotations Panel. The measurement is deleted.

To delete a distance measurement in the L-Mode:

1. Click on a distance measurement in the L-Mode to select it. A white box appears around the label.
2. Press the <Del> key on the keyboard.

The distance measurement is deleted.

To delete a text callout:

1. Click on a text callout to select it. A white box appears around the text callout.
2. Press the <Del> key on the keyboard.

The text callout is deleted.

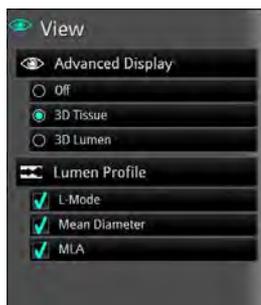
## Deleting All Measurements and Text Callouts



Click to delete all measurements, calculations, and text callouts from this recording or still image.

**NOTE:** Automatically generated **Lumen Contours** area measurements are not deleted. If you have used the **Lumen Profile** function to calculate the minimum lumen area, that measurement is not deleted.

## Lumen Profile Display Option



The system automatically creates a trace of the lumen contour on each frame.

**NOTE:** Lumen Profile is not available with still images or stationary recordings.

L-Mode view must be turned on in order to use the Lumen Profile display option.

Remove the check from the Measurements checkbox to hide the lumen contour trace and corresponding measurements.

The Lumen Profile display opens with Minimum Lumen Area (**MLA**) and a percent stenosis graph (**%AS** or **%DS**) turned on:


- **MLA** displays the Minimum Lumen Area controls and values for this recording.
- **%DS** shows how the mean diameter changes along the L-Mode. It can also display the **%DS** at **MLA** and proximal and distal reference points.
- **%AS** shows how the lumen area changes along the L-Mode. It can also display the **%AS** at **MLA** and proximal and distal reference points.

**NOTE:** For information on displaying extended MLA info, see [Table 10-12 on page 10-21](#).

## Lumen Profile Display With MLA Controls Overview

**CAUTION:** It is the user’s responsibility to confirm the lumen contour on each frame, and to make adjustments if necessary. Red frames indicate low confidence in the detected contour. If the MLA frame is in the vicinity of a low confidence region, the system displays “?.??” for the MLA values. If the MLA values are “?.??,” the contour on the MLA frame must be reviewed, edited if necessary, and accepted before it can be displayed.

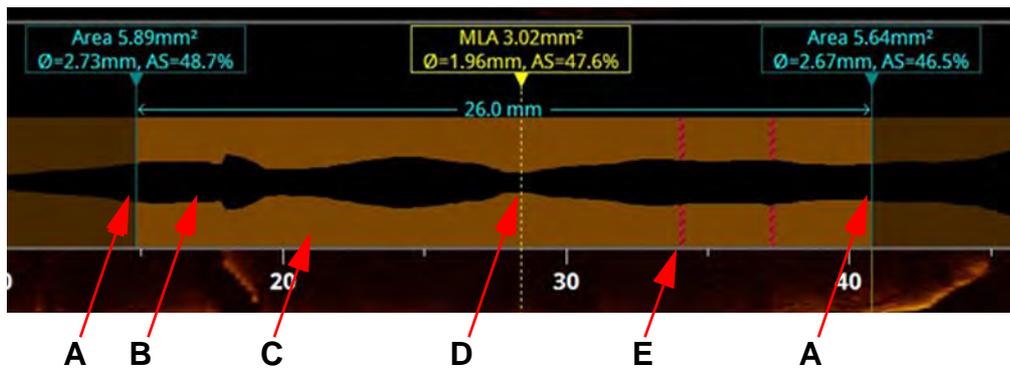Table 7-2: MLA Controls



**A Distal and Proximal Reference Frames :** Move these to set the range for MLA calculation. The system’s search for a minimum lumen area occurs on frames between the distal and proximal reference frames.

See [Table 7-3 on page 7-22](#) for more information on calculations. See [Table 10-12 on page 10-21](#) for showing extended MLA info.

- If %DS is selected in the Lumen Profile, the reference frames and the MLA frame display the percent diameter stenosis at each end of the range.
- If %AS is selected, the reference frames and the MLA frame display the percent area stenosis at each end of the range.

**B** The lumen area is colored black.

**C** Where the system has high confidence in the contour of the lumen area, or where the contour has been accepted by the user, the section is colored brown.

**D Calculated MLA :** The dashed line indicates the position of the minimum lumen area between the distal and proximal reference frames (A).

**E** Where the system has low confidence in the contour of the lumen area, the section is colored red. These frames are not considered in the MLA search. If these frames are within the range where MLA is calculated, you must go to the MLA frame and confirm the contour.

**Measurements and Annotations**  
**Lumen Profile Display Option**


---

Table 7-3: Lumen Profile %AS and %DS Calculations

	<b>When Lumen Profile set to %DS...</b>	<b>When Lumen Profile set to %AS...</b>
Stenosis calculations of the MLA frame compared to average of distal and proximal references (value listed above the MLA frame):	$\%DS = 100 * (D_{REF} - D_{MLA}) / D_{REF}$ <p>where <math>D_{REF}</math> is the average of the Distal Reference's mean diameter and the Proximal Reference's mean diameter. <math>D_{MLA}</math> is the MLA's mean diameter.</p>	$\%AS = 100 * (A_{REF} - A_{MLA}) / A_{REF}$ <p>where <math>A_{REF}</math> is the average of the Distal Reference's area and the Proximal Reference's area. <math>A_{MLA}</math> is the Minimum Lumen Area.</p>
Stenosis calculations of the MLA frame compared to individual distal and proximal references (values listed above the D or P reference frames):	$\%DS = 100 * (D_{D \text{ or } P} - D_{MLA}) / D_{D \text{ or } P}$ <p>where <math>D_{D \text{ or } P}</math> is the Distal Reference's mean diameter or the Proximal Reference's mean diameter. <math>D_{MLA}</math> is the MLA's mean diameter.</p>	$\%AS = 100 * (A_{D \text{ or } P} - A_{MLA}) / A_{D \text{ or } P}$ <p>where <math>A_{D \text{ or } P}</math> is the Distal Reference's area or the Proximal Reference's area. <math>A_{MLA}</math> is the Minimum Lumen Area.</p>

## 3D Display Option

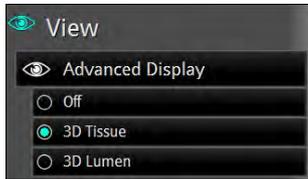Click the **Advanced Display** button to create a 3D representation of the recording.



Clicking the **Advanced Display** button toggles-through the options: **Off**, **3D Tissue**, and **3D Lumen**. You can also select the desired option by clicking the button next to it.

**NOTE:** You cannot add or change measurements or annotations in the cross-section view or 3D Display while 3D Display is on. You can add annotations and distance measurements to the L-Mode view only. 3D Display is not available with still images or stationary recordings.

## Measurements and Annotations


### 3D Display Option

---

### 3D Tissue Controls

Table 7-4: 3D Tissue Controls

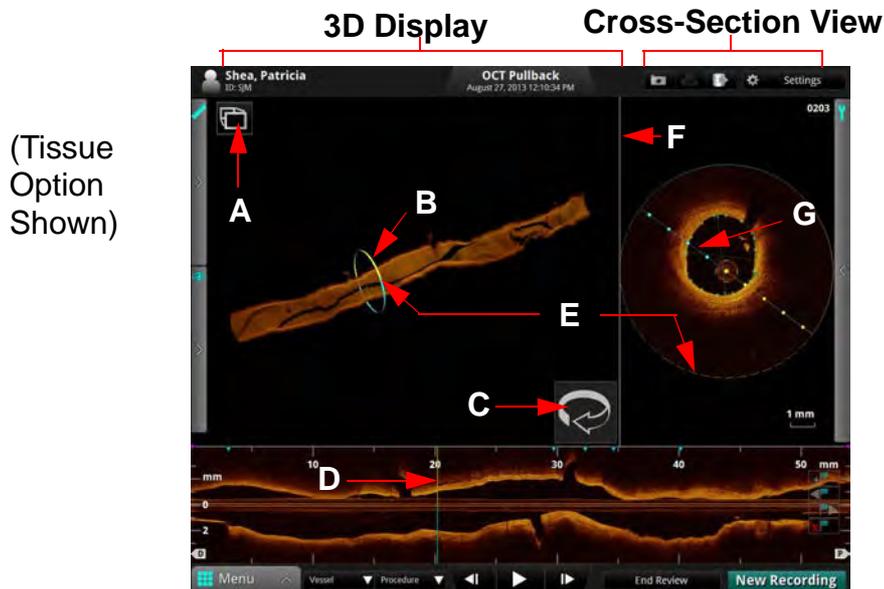---



- A View Mode** : Click to toggle between **Window** mode (shown above) and **Full Screen Mode**.
  - B Current Frame Indicator (3D View)** : Scroll the mouse wheel (or twist the **Navigation Controller**) to change the cut plane as shown in the L-Mode.
  - C Cut Plane Rotation Hotspot** : Place the cursor over the Hotspot and rotate the mouse wheel to change the cut plane as shown in the L-Mode.
  - D Current Frame Indicator (L-Mode View)** : Click and drag to change the frame shown.
  - E** The solid lines in the cross-sectional view represent the rendered half of the vessel image. The dashed lines in the cross-sectional view represent the open, or un-rendered, half of the vessel image. The blue and yellow colors are for location referencing among the views.
  - F** Click and drag the divider bar side to side to change the size of 3D Display versus cross-section view.
  - G Cut Plane Indicator** : The cut-plane is shown as a solid line in the cross-sectional view Click and drag this to change the lateral view shown in the L-Mode display.
- NOTE:** To zoom you can either left click in the 3D Display window for Quick Zoom factors or right-click and drag up or down, left or right to adjust zoom.
-

### 3D Display with Lumen

**NOTE:** If the 3D Lumen setting is turned on, the 3D Tissue setting is turned off, and vice-versa. Selecting one deselects the other.

The **3D Lumen** setting adds a 3D representation of the lumen contours drawn on each frame.



Figure 7-13: 3D Display with Lumen

**NOTE:** If you have the **MLA** setting turned on under the **Lumen Profile** display and you have selected **3D Lumen** (shown) or **3D Tissue**, the distal, proximal, and MLA frames are displayed.

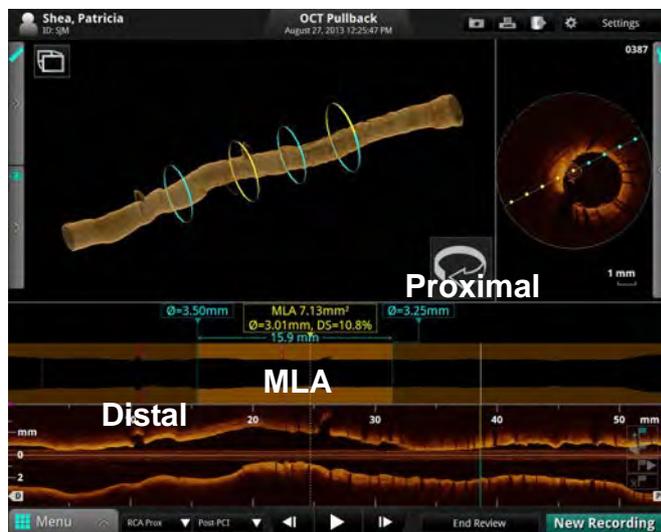


Figure 7-14: MLA Frames in 3D

### **Limitations of 3D Display**

Please be aware of the following limitations of the 3D Display Option:

- Although a vessel may curve, the 3D Display view always appears straight because of the limitations of reconstruction.
- Artifacts caused by the relative motion of the catheter and the vessel often result in a saw-toothed appearance of the reconstruction and can lead to misinterpretation by inexperienced users.
- Shortening or lengthening artifacts in the 3D Display reconstruction may occur due to the relative motion of the imaging catheter with respect to the coronary artery caused by the patient's heart motion.

# Exporting, Importing, and Managing Files

---

8

**CAUTION:** Please note St. Jude Medical makes no representation or warranty that use of the OPTIS INTEGRATED System complies with applicable privacy, security and confidentiality laws, but encourages you to assess your own risk as you use, disclose, control, process or transfer patient health information with the OPTIS INTEGRATED System.

## **Compatible Transfer Media and USB Devices**

The following sections list the supported media formats.

If you experience problems with a specific type of CD/DVD or USB device, contact St. Jude Medical for recommendations.

### **Optical Media**

The system can export files through the CD/DVD drive. [Table 8-1](#) lists the supported disc formats and describes whether they can be erased and/or appended.

Table 8-1: Optical Media Characteristics

<b>Media</b>	<b>Capacity</b>	<b>Can Erase</b>	<b>Can Append</b>
CD-R	737 MB		X
CD-RW	737 MB	X	X
DVD+R	4.7 GB		X
DVD+RW	4.7 GB	X	
DVD+R Dual Layered	8.5 GB		X
DVD-R	4.7 GB		X
DVD-RW	4.7 GB	X	
DVD-R Dual Layered	8.5 GB		
DVD-RAM	4.7 GB	X	

If you try to export files to a CD-RW that already has data on it, an alert message appears. Click **Yes** to export your files to this disc, or click **No** to cancel the export.

If you try to export files to a DVD+RW, DVD-RW, or DVD-RAM that already has data on it, an alert message appears. Click **OK** to erase the disc and export your files, or click **Cancel** to cancel the export.

If you try to export files to a disc that cannot be formatted or appended (such as DVD+R DL and DVD-R DL), an alert message appears. Insert a new disc to continue, or click **Cancel** to cancel the export.

### **USB Connected Media**

The system can export files through the USB port. Any USB hard drive or Flash Drive that meets the USB 2.0 Specification and is supported by Windows 7 SP1 may be used.

**WARNING:** Inside the catheterization lab only port-powered USB drives may be connected to the USB port. Connecting externally powered devices to

**the USB port in the patient vicinity may compromise electrical isolation and cause patient injury.**

**NOTE:** Outside the catheterization lab, IEC 60950-compliant, externally powered USB hard drives may be connected to the USB port.

## File Formats

Files can be exported in native file format (raw OCT format), a standard graphic file format (standard format), or DICOM format. Exported files can be saved on a CD/ DVD or external USB device, or exported to a DICOM storage server. You can choose whether to delete files after transfer or keep them on the system.

### About Native (Raw OCT) Format

If you export in native (raw OCT) file format, every feature of the OCT file will be exported, and the files can be imported into another ILUMIEN or OPTIS INTEGRATED System or an Offline Review Workstation (ORW) and reviewed and manipulated there. An exported OCT file contains exactly the same data as the original file, including any measurements and annotation additions, all patient information associated with each file, and system diagnostic information to help diagnose possible image quality problems.

Note the following:

- Native (raw OCT) files can be reviewed and manipulated only with an ILUMIEN System, an OPTIS INTEGRATED System, or an ORW. They cannot be accessed with any other systems or software.
- Native (raw OCT) files are large; exporting to CD/DVD may require many discs.

### About DICOM Format

For a multiple-frame recording using DICOM format stored on a network server, the maximum resolution is 800 x 800.

## About Standard Format

If you export in standard file formats (AVI, compressed AVI, or Multi-page TIFF for recordings; JPEG, TIFF, or BMP for images), the images can be used in computer applications outside the system, but cannot be imported into an ILUMIEN System, an OPTIS INTEGRATED System, OPTIS MOBILE System, or an ORW.

**NOTE:** For standard format multiple-file exports, the same limitations apply as described in [“Image Format and Size in Standard Formats”](#). Just as when exporting a single file, you can choose the size and format.

## **Image Format and Size in Standard Formats**

When exporting files in standard formats, the system gives you several choices of format and file resolution (size). When making these choices, keep in mind the following:

### **File Size**

- File size is dependent on resolution; the lower the resolution, the smaller the file.
  - Decreasing resolution makes the resulting image grainier, but no quantitative data, annotations or measurements are lost.
  - Increasing resolution results in higher image quality in the exported file. For a file that is a single frame (from a captured image, a paused image, or from a recording that has been edited down to only one frame), the maximum resolution for export is 2048 x 2048.
  - For a multiple-frame recording, the maximum resolution is 1024 x 1024.
- It is best to select the resolution that will actually be used in the intended application. When images are enlarged or reduced by the system, the maximum useful information is preserved. External applications may not apply the same diligence when enlarging or reducing an image.

### **Standard File Format**

- The system can export recordings as either AVI, compressed AVI, or Multi-page TIFF format.
  - Most computers include standard players that will play AVI format files.
  - The Compressed AVI and Multi-page TIFF formats require that special players for these formats be installed on the computer you will use to review the files.
- The system can export still images as either JPEG, TIFF, or BMP.
  - TIFF and BMP files are high resolution bit-map files with large file size.
  - JPEG files are compressed image files designed especially for viewing in a computer application. They are significantly smaller than TIFF or BMP files, and are typically the best choice for slide shows or other applications where the files will be viewed on a computer monitor.

## Exporting Files During a Review

During review, you can use the **Export** button to export the current recording or still image to Native (Raw) format, DICOM format, or standard file formats.

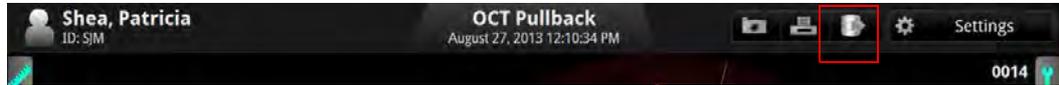


Figure 8-1: Export Button (OCT)

**NOTE:** The **Export** button is unavailable while a recording is playing. To export from a recording, pause the recording before clicking the **Export** button.

All edits, including measurements, calculations, annotations and zoom, are included with the recording or still image when it is exported (see [Chapter 7 “Measurements and Annotations”](#)).

**CAUTION:** If you are exporting to standard formats, make all measurements **BEFORE** exporting. Using non-OCT software to measure standard format images will not produce accurate measurements.

**CAUTION:** Do not use images that have been exported to JPEG or Compressed AVI formats for clinical decision making. These formats use compression methods that may degrade the image quality.

**NOTE:** Editing images stored on the system does not change the underlying image data in any way. All unedited data taken during the procedure is always preserved and can be reviewed without the changes.

**NOTE:** If the system displays an error message during export, a restart may be required. If necessary, click **OK** to restart the system. If the error condition persists, export the recording from the Patient Summary menu, or reduce the resolution for the export of the recording. See [“Exporting Files from the Patient Summary Menu”](#) on page 8-14 for more information.

#### Exporting Files in Native (Raw) Format

1. Click the **Export** button.

The **Export Wizard** opens (see [Figure 8-2](#)).



Figure 8-2: The Export Wizard - Step 1

2. Click to select **Native (Raw)**, and click **Next**.
3. If you need to remove patient identifying information from the exported files, add a check to the **Anonymize** checkbox.
  - If you have added a check to the **Anonymize** checkbox, you can also add an alternate means of identifying the patient. Click the **Alternate PIDs...** button to open the **Define Alternate Patient ID** menu:
    - Click in the **Alternate Patient ID** column and type an alternate ID tag for this patient.



Figure 8-3: Define Alternate Patient ID Menu

- Click the **Auto-Generate** button to add an alternate ID generated by the system. Click the **Reset** button to clear the ID generated by the system.

- Click **OK** to approve the change and close the menu, or **Cancel** to return to previous settings and close the **Define Alternate Patient ID** menu.
4. Click to set the **System File Options** for the file(s) that you are exporting:
- **Leave Unchanged** - the files are exported, but the original files in the System are not changed.
  - **Mark as Archived** - the files are exported, and marked as archived.
  - **Remove when Complete** - the files are deleted from the System once the export is complete.
5. Click **Next** to approve the settings and open the next page of the **Export Wizard**.

**NOTE:** Click **Back** to return to the previous page, or click **Cancel** to exit the Export Wizard without exporting files.

6. Click to select an **Export Destination**:

**NOTE:** If an output device is not available, that option cannot be selected.

**NOTE:** Check that the **Free** space in the Drive Capacity section is large enough to save the files that are being exported.

- Click **CD/DVD** to export to a CD/DVD.

If the inserted CD/DVD is blank, you can edit the name in the **Volume Label** field.

- Click **External Drive** to export to an external USB device.

If necessary, click the **External Drive** selection box to select the correct drive location from the drop-down menu.

Click the **Ellipsis (...)** button to the right of the **Path** selection box to browse for a specific folder on the device.

7. Click **Next** to approve the settings and open the next page of the **Export Wizard**.

The next page of the **Export Wizard** opens. The menu shows a summary of Export Information.

8. Click **Export** to export the selected files.

The export begins. If necessary, click **Cancel** to stop the export in progress.

## Exporting, Importing, and Managing Files

### Exporting Files During a Review

---

9. When complete, click **Done** to close the **Export Wizard**.

### Exporting Files in DICOM Format

1. Click the **Export** button.

The **Export Wizard** opens (see [Figure 8-2](#)).

2. Click to select **DICOM**, and click **Next**.

**NOTE:** If you are reviewing a single frame instead of a recording, skip to [Step 5](#)

3. Click to select the material to be exported:

- **Pullback**
- **Current Frame**
- **Bookmarked Frames** - This choice is unavailable if there are no bookmarked frames in the recording.

4. Click **Next** to approve the settings and open the next page of the **Export Wizard**.

**NOTE:** Click **Back** to return to the previous page, or click **Cancel** to exit the Export Wizard without exporting files.

5. If you need to remove patient identifying information from the exported files, add a check to the **Anonymize** checkbox.

6. Click to set the **System File Options** for the file(s) that you are exporting:

- **Leave Unchanged** - the files are exported, but the original files in the System are not changed.
- **Mark as Archived** - the files are exported, and marked as archived.

7. Click **Next** to approve the settings and open the next page of the **Export Wizard**.

8. Click to select an **Export Destination**:

**NOTE:** If an output device is not available, that option cannot be selected.

**NOTE:** Check that the **Free** space in the Drive Capacity section is large enough to save the files that are being exported.

- Click **CD/DVD** to export to a CD/DVD.

If the inserted CD/DVD is blank, you can edit the name in the **Volume Label** field.

- Click **External Drive** to export to an external USB device.

If necessary, click the **External Drive** selection box to select the correct drive location from the drop-down menu.

Click the **Ellipsis (...)** button to the right of the **Path** selection box to browse for a specific folder on the device.

- Click **Remote Store** to export to a network DICOM Server.

**NOTE:** The **Remote Store** option is enabled only if the DICOM network connection has been configured on the **DICOM** tab of the **Setup** dialog box (see [Chapter 10 “User Interface Reference”](#) for more information on DICOM settings).

9. Click **Next** to approve the settings and open the next page of the **Export Wizard**.

The next page of the **Export Wizard** opens, showing a summary of Export Information.

10. Click **Export** to export the selected files.

The export begins. If necessary, click **Cancel** to stop the export in progress.

11. When complete, click **Done** to close the **Export Wizard**.

#### Exporting Files in Standard Formats

**NOTE:** 3D reconstructions cannot be exported in Standard formats.

1. Click the **Export** button.

The **Export Wizard** opens (see [Figure 8-2](#)).

2. Click to select **Standard**, and click **Next**.

**NOTE:** If you are reviewing a single frame instead of a recording, skip to [Step 5](#).

3. Click to select the material to be exported:

- **Pullback**
- **Current Frame**
- **Bookmarked Frames** - This choice is unavailable if there are no bookmarked frames in the recording.

4. Click **Next** to approve the settings and open the next page of the **Export Wizard**.

**NOTE:** Click **Back** to return to the previous page, or click **Cancel** to exit the Export Wizard without exporting files.

5. Click to set the resolution and file format.

6. Click **Next** to approve the settings and open the next page of the **Export Wizard**.

7. Verify the names of the files that are being exported. If necessary, double-click a filename and edit it using the keyboard.

8. Click **Next** to approve the settings and open the next page of the **Export Wizard**.

9. Click to select an **Export Destination**:

**NOTE:** If an output device is not available, that option cannot be selected.

**NOTE:** Check that the **Free** space in the Drive Capacity section is large enough to save the files that are being exported.

- Click **CD/DVD** to export to a CD/DVD.

If the inserted CD/DVD is blank, you can edit the name in the **Volume Label** field.

- Click **External Drive** to export to an external USB device.

If necessary, click the **External Drive** selection box to select the correct drive location from the drop-down menu.

Click the **Ellipsis (...)** button to the right of the **Path** selection box to browse for a specific folder on the device.

10. Click **Next** to approve the settings and open the next page of the **Export Wizard**.

The next page of the **Export Wizard** opens. The menu shows a summary of Export Information.

11. Click **Export** to export the selected files.

The export begins. If necessary, click **Cancel** to stop the export in progress.

12. When complete, click **Done** to close the **Export Wizard**.

## Exporting Files from the Patient Summary Menu

You can export files from the **Patient Summary** menu.

1. In the **Patient Summary** menu, add a check to the checkbox of each file that you want to export.

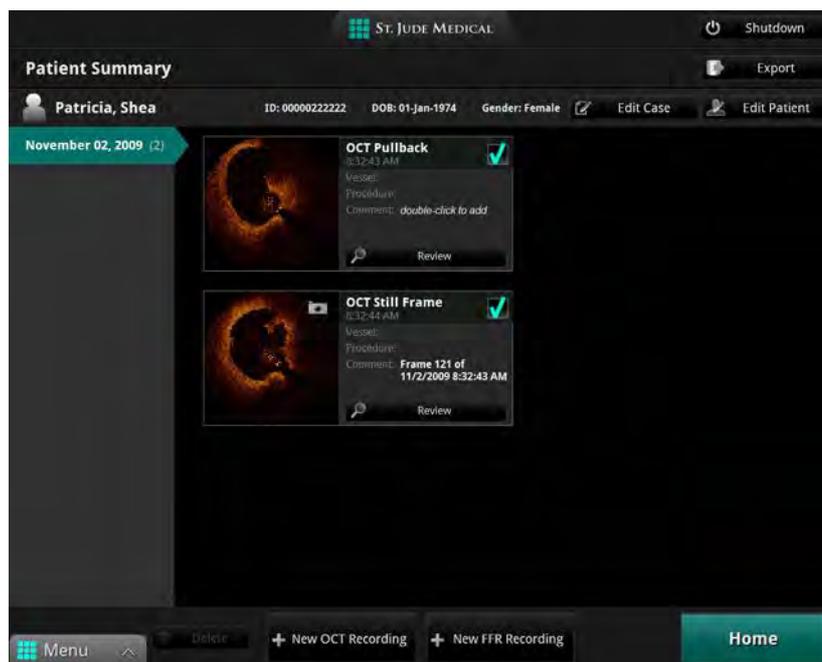


Figure 8-4: Highlighted Records

A white border appears around each record that is selected.

**NOTE:** In the **Patient Summary** menu, you can select multiple files from the same date, but not from different dates. Within the **Export Wizard**, you can select additional files from different dates and different patients.

**NOTE:** To deselect a file, clear the check from the checkbox. The check and white border disappear to indicate that the file is no longer selected.

2. Click the **Export** button at the top of the screen.

The **Export Wizard** opens (see [Figure 8-2](#)).

3. Click to select an export format, and click **Next**.

The next page of the **Export Wizard** opens.

4. To filter the list, click the **All**, **Unarchived**, or **New since last archive** button.
5. If necessary, add a check to the checkbox of any other recording to be exported.
  - Add a check to the checkbox of patient to select all recordings for that patient.
  - Add a check to the checkbox of the first entry in the list to select all recordings.
6. Click **Next** to approve the list of files for export and open the next page of the **Export Wizard**.
  - If you are exporting in Native (Raw) format, refer to [Step 3](#) in “[Exporting Files in Native \(Raw\) Format](#)” on page 8-8 to continue the export.
  - If you are exporting in DICOM format, refer to [Step 5](#) in “[Exporting Files in DICOM Format](#)” on page 8-10 to continue the export.
  - If you are exporting in standard formats, refer to [Step 5](#) in “[Exporting Files in Standard Formats](#)” on page 8-12 to continue the export.

## Using Exported Standard Format Recordings

The system can export recordings as either AVI, compressed AVI, or Multi-page TIFF format. Most computers include standard players that will play AVI format files. However, the Compressed AVI format and the Multi-page TIFF formats require that you have special players installed on your computer to review the files.

### St. Jude Medical DICOM Viewer

When an image file is exported to a CD/DVD in the DICOM format, the St. Jude Medical DICOM Viewer (shown in Figure 8-5) is also included on the CD/DVD. This DICOM CD/DVD can then be freely distributed and used with any Microsoft Windows computer (the CD/ DVD is configured to automatically launch the viewer when it is inserted into the drive). The St. Jude Medical DICOM Viewer may only be used to view St. Jude Medical generated DICOM images on a St. Jude Medical DICOM CD/DVD. If this CD/DVD is to be used in a public forum, the Anonymous option should be used when exporting the image using the system.



Figure 8-5: St. Jude Medical DICOM Viewer - Image View

## Exporting, Importing, and Managing Files Using Exported Standard Format Recordings

The St. Jude Medical DICOM Viewer window contains three different sections:

- **Viewer Controls** - Provides control over the viewer and the active image.
- **Directory Tree** - This tree lists the patient, study, series, and image hierarchy contained in the DICOMDIR file on the CD/DVD. Individual images may be loaded by clicking on the image item in the tree.
- **File Viewer** - Provides two tabs:
  - **Image** tab - Displays the image currently selected in the DICOM Directory Tree (see Figure 8-5). If the image is a recording, playback may be controlled using the playback buttons in the Viewer Controls at the top of the viewer. The image may also be zoomed in or out using the **Zoom** feature in the toolbar.
  - **Attributes** tab - Displays a list of all DICOM modules and attributes contained in the currently selected image (see Figure 8-6).

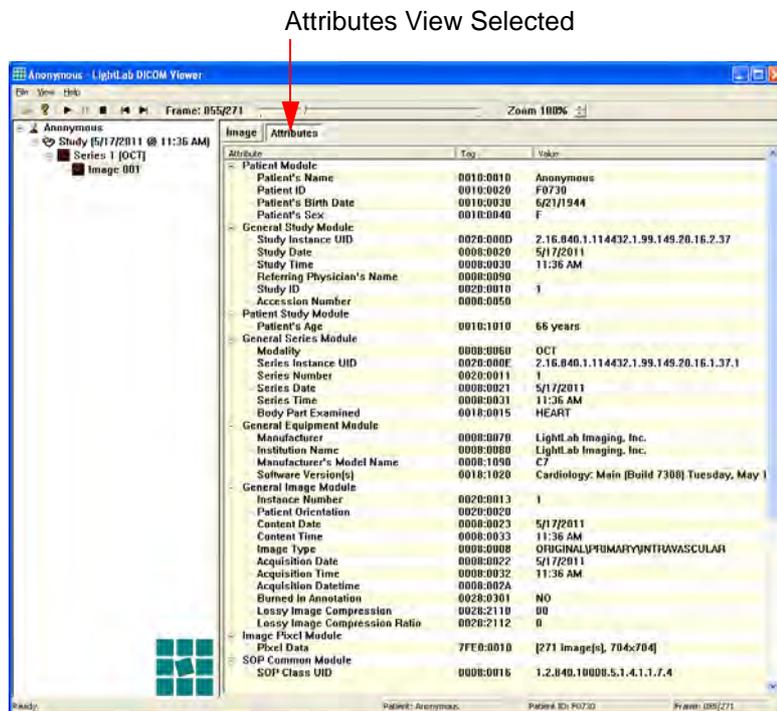


Figure 8-6: St. Jude Medical DICOM Viewer - Attributes View

## Exporting, Importing, and Managing Files Using Exported Standard Format Recordings

The St. Jude Medical DICOM files are compliant with the DICOM Standard (PS 3-2008) and use the Secondary Capture Multi-Frame Image Storage IOD (Information Object Definition). The actual attributes contained in the file are listed in [Table 8-2](#).

Table 8-2: DICOM File Attributes

Module	Attribute
Patient	Patient's Name
	Patient ID
	Patient's Birth Date
	Patient's Gender
General Study	Study Instance UID
	Study Date
	Study Time
	Referring Physician's Name
	Study ID
	Accession Number
	Patient Study
General Series	Modality
	Series Instance UID
	Series Number
	Series Date
	Series Time
	Body Part Examined
	General Equipment
Institution Name	
Manufacturer Model Name	
Software Version(s)	
General Image	Instance Number
	Patient Orientation
	Content Date
	Content Time
	Image Type
	Acquisition Date
	Acquisition Time
	Acquisition DateTime
	Burned In Annotation
	Image Comments

## Exporting, Importing, and Managing Files Using Exported Standard Format Recordings

Table 8-2: DICOM File Attributes (*continued*)

<b>Module</b>	<b>Attribute</b>
Image Pixel SOP Common	Lossy Image Compression
	Lossy Image Compression Ratio
	Pixel Data
	SOP Class UID
	SOP Instance UID
	Instance Creation Date
	Instance Creation Time
	Timezone Offset from UTC
	Instance Number
	Frame Increment Pointer
Multi-frame <sup>1</sup>	Sequence of Regions
Region Calibration	Frame Time
Cine <sup>1</sup>	Conversion Type
SC Equipment	Modality
	SC Device Manufacturer
	SC Device Model Name
	SC Device Software Version
SC Multi-frame Image <sup>1</sup>	Burned In Annotation
	Frame Increment Pointer

1. Multi-frame files only.

## Importing Files from a CD/DVD or USB

1. Connect an external database to the ILUMIEN OPTIS System:
  - Insert a CD/DVD containing a database into the system.
  - Connect an external hard drive to the system.
2. In the **Select Patient** menu, click the Import button.

The **Select Database** menu opens.

3. Select the database to import, and click OK.
4. The **Import Database** menu opens (Figure 8-7).

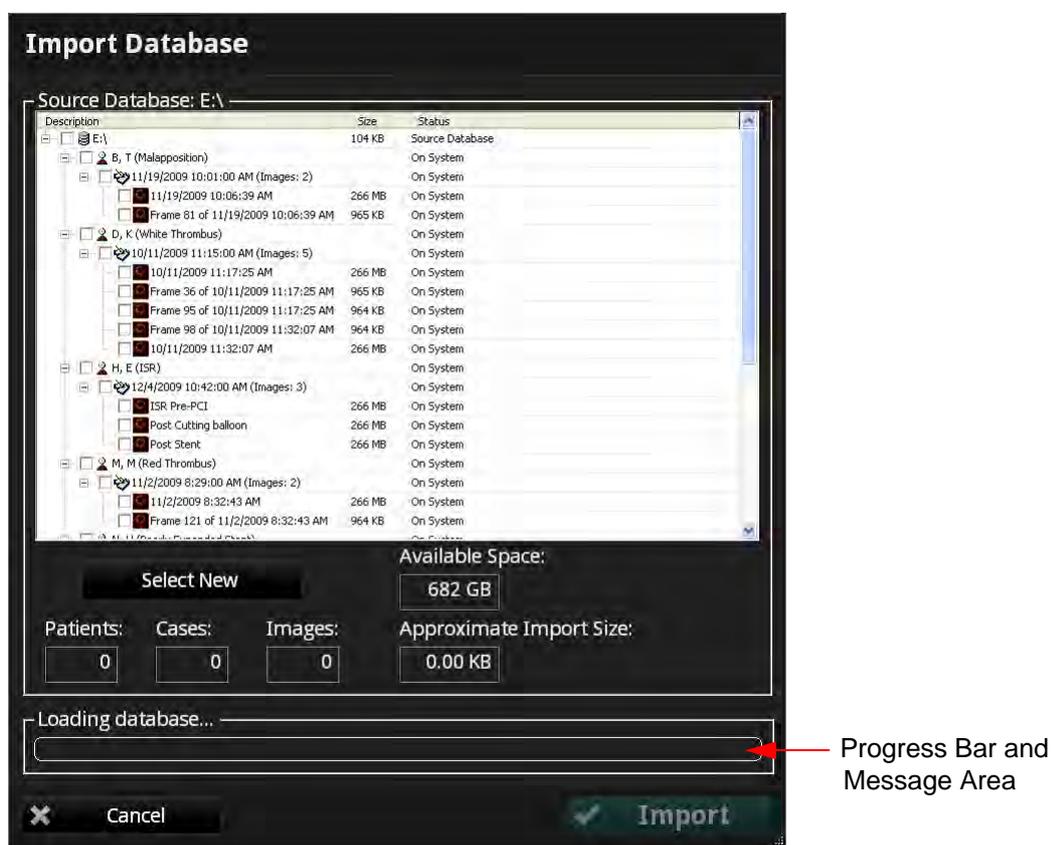


Figure 8-7: Import Database Menu

Table 8-3: Import Database Menu Options

<b>Source Database</b>	<p>List of files which will be imported. Each file is displayed in one row, and the patient name, patient ID, image creation date &amp; time, comment, and status are displayed in columns. The list is initially sorted by ascending image creation date &amp; time and can be sorted by any of the columns in ascending or descending order by left-clicking once or twice on the desired column header.</p> <p>As each file is imported, its <b>Status</b> message is updated:</p> <ul style="list-style-type: none"><li>• <b>Importing</b> - file is currently being imported.</li><li>• <b>On System</b> - file has been copied to the system.</li><li>• <b>On Archive Media</b> - file is on the current archive media and has not yet been imported.</li><li>• <b>Absent</b> - file is not located on the current archive media</li><li>• <b>Failed</b> - attempt to import the file failed.</li></ul>
<b>Select New</b>	<p>Click to select only recordings that are not present in the ILUMIEN OPTIS System.</p>
Progress Bar and Message Area	<p>The Progress Bar indicates the status of the file transfer. The Message Area displays information about the file and the transfer.</p>
<b>Import</b>	<p>Import the selected files.</p>
<b>Cancel</b>	<p>Cancel the import.</p>

---

An import may require installation of several CD/DVDs. When a new disc is required, the system ejects the current disc and displays a **New Disc Request** alert.

## Deleting Files

You can delete files by any of the following methods:

- Exporting images with the **Remove When Complete** function turned on (see [“Exporting Files in Native \(Raw\) Format”](#) on page 8-8).
- Select and delete files from the **Patient Summary** menu.
- Select and delete files from the **Database** menu.

**CAUTION:** Once files are deleted, they cannot be restored. After files have been deleted, they can only be imported back to your system from your archived copies.

### Deleting Files from the Patient Summary Menu

1. In the **Patient Summary** menu, add a check to the checkbox of each file that you want to delete (see [Figure 8-4 on page 8-14](#)).

A white border appears around each record that is selected.

**NOTE:** In the **Patient Summary** menu, you can select multiple files from the same date, but not from different dates.

**NOTE:** To deselect a file, clear the check from the checkbox. The check and white border disappear to indicate that the file is no longer selected.

2. Click the **Delete** button at the bottom of the screen.

A prompt asks you to confirm that you want to delete the selected files.

3. Click **Yes** to delete the files, or click **No** to cancel the deletion and return to the **Patient Summary** menu.

## Deleting Files from the Database Menu

1. Click the **Menu** button and select **Setup**. The **Setup** dialog box opens.
2. Click the **Database** button to open the **Database** menu (see [Figure on page 10-7](#)).
3. Add a check to the checkbox of any recording to be deleted.
  - Add a check to the checkbox of patient to select all recordings for that patient.
  - Add a check to the checkbox of the first entry in the list to select all recordings.
4. Click the **Delete** button.

A prompt asks you to confirm that you want to delete the selected files.

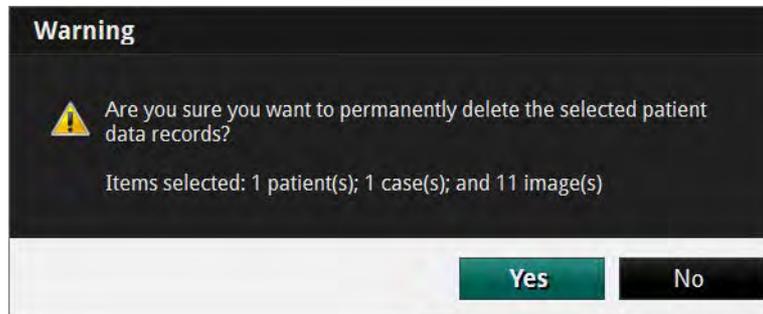


Figure 8-8: Deletion Warning Alert

5. Click **Yes** to delete the files, or click **No** to cancel the deletion and return to the **Database** menu.

## Transfer and Import Messages

The system displays a number of different messages to inform you of problems that may arise during transfer and import operations.

Table 8-4: Transfer Messages

<b>Error Message</b>	<b>Cause</b>	<b>Action</b>
The CD/DVD drive is empty. Please insert a writable CD/DVD, or click Cancel to stop the transfer.	No CD/DVD is in the drive, though you have attempted to transfer files.	To continue with the export, insert a blank CD/ DVD or a previously used CD/DVD on which there is sufficient space for the transfer. To stop the transfer, click <b>Cancel</b> .
None of the selected files will fit on the disc. Please insert a different CD/DVD, or click Cancel to stop the transfer.	Insufficient space is available on the destination CD/DVD and NONE of the selected files can be transferred to this disc.	Insert a new writable CD/DVD into the drive and close the drive to transfer or click <b>Cancel</b> to stop the transfer and file deletion.
Some files did not fit on the disc. Please insert another writable CD/DVD to continue the transfer, or click Cancel to stop the transfer.	Insufficient space is available on the destination CD/ DVD for all of the selected files, but SOME of the files have been transferred.	To continue the transfer on the new disc, insert another writable CD/DVD into the drive and close the drive to continue the transfer. To interrupt the transfer, click <b>Cancel</b> . Some files may be on the first disc, some on the second. If you select Cancel, the CD/ DVD will contain the files previously transferred, but the remainder will not be transferred and the uncopied files will remain on the hard drive.
Cancelling the transfer may damage the disc and render previously written files unreadable. Are you sure you want to cancel the transfer?	You have attempted to interrupt a transfer in process by clicking <b>Cancel</b> .	Click <b>No</b> to continue the transfer. Click <b>Yes</b> to stop the transfer and risk losing the files already transferred.
An invalid volume name was specified.	You have typed a character that cannot be used in a volume name.	Retype the name with acceptable characters.
Transfer failed, the destination disk is full.	The USB drive is full.	Connect a different USB drive.

Table 8-4: Transfer Messages (*continued*)

<b>Error Message</b>	<b>Cause</b>	<b>Action</b>
Transfer failed, the destination disc is not blank.	The CD/DVD already contains files and cannot be appended.	Insert a new writable CD/DVD into the drive and close the drive to transfer or click <b>Cancel</b> to stop the transfer.
Transfer failed, the source or destination disk was removed.	The USB drive is disconnected.	Reconnect the USB drive.
Please insert the next CD to continue importing the database, or click Cancel to stop importing.	The current disk is full and has been ejected.	Insert another CD/DVD in the series or click <b>Cancel</b> to stop the import.
Duplicate CD. Please insert the next CD to continue importing the database, or click Cancel to stop importing.	You have inserted a CD/DVD containing files that have already been imported.	Insert another CD/DVD in the series. or click <b>Cancel</b> to stop the import.
Database not found. Please insert the next CD to continue importing the database, or click Cancel to stop importing.	The currently inserted CD/DVD does not include a valid OCT database file.	Insert a CD/DVD which contains a valid OCT database file (OCT.mdb or OCT.dbf) or click <b>Cancel</b> to stop the import.
Incorrect database. Please insert the next CD to continue importing the database, or click Cancel to stop importing.	The currently inserted CD/DVD includes a valid OCT database file, but is not in the same series as the previously inserted CD/DVD(s).	Insert a CD/DVD from the same series as the last CD/DVD or click <b>Cancel</b> to stop the import.
Database import not complete.	You have clicked <b>Cancel</b> , so the Import operation is not complete.	Click <b>OK</b> to continue system operations. Some of the files may have been partially, but not fully copied into the system. If you try to open one of these files, the system will display an error message.
Please insert the next CD to continue importing the database, or click Cancel to stop importing.	Additional files exists on another disc.	Insert the next disc. The message is automatically closed and import continues.

## Duplicate File Name Messages

Duplicate file name messages warn you when a file has been selected for export with the same name as one already saved. The system does not allow any files to be overwritten.

**NOTE:** Do not save an image using the same name used for a previous file. St. Jude Medical recommends that you label media containing exported files with descriptive information, including the procedure or export date. To resolve these errors when copying to CD/DVDs, insert another CD/DVD into the drive and close the drive. File transfer resumes automatically. If the transfer is cancelled, the uncopied files will remain on the hard drive.

Table 8-5: Duplicate File Name Messages

Message	Cause	Action
Duplicate files could not be written. Please insert a different CD/DVD to continue the transfer, or click Cancel to stop the transfer	SOME of the file names in an export set are already found on the installed CD/DVD.	Files with non-identical names will be written to the current CD/DVD; files with identical names will not be written. To save the files with the duplicate file names, insert a new CD/DVD and continue with the transfer. To stop the transfer, click <b>Cancel</b> .
The disc already contains files with the same names. Please insert a different CD/DVD to continue the transfer, or click Cancel to stop the transfer.	ALL the file names in an export set are already found on the installed CD/DVD.	To save these files with these file names, insert a new CD/DVD and continue with the transfer. To stop the transfer, click <b>Cancel</b> . This message may be displayed when the <b>Delete after transfer</b> checkbox in the <b>Manage Exported Files</b> menu is unchecked, leading to re-exporting the same files.

## Database Statistics

You can check the size and statistics of your system's database using the **Database Maintenance** button in the **Setup** dialog box. See "[Setup - Database/Maintenance Menu](#)" on [page 10-9](#) for more information.

## Importing Patient Information From a DICOM Worklist or Storage Server

### From DICOM Worklist

To import patient information From a DICOM Worklist, complete the following steps:

1. From the Select Patient menu, click the **Add Patient** button. The **Add Patient - Step 1** guidance displays.



Figure 8-9: Add Patient - Step 1 (Worklist)

2. Select **DICOM Worklist Items**. Click the **Next** button. The **Add Patient - Step 2** guidance displays.

**+ Add Patient - Step 2 of 3**  
Specify modality workload query information.

Modality Worklist Query Information Server: Siemens\_DMWL Test Connection

Patient Name	Requested Procedure ID
Patient ID	Modality: IVQCT
Accession Number	<input checked="" type="checkbox"/> Scheduled Procedure Step Start Date
Scheduled Station AE	From: 01-Nov-2013 To: 11-Feb-2014

Cancel Back Next

Figure 8-10: Add Patient - Step 2

3. Enter a valid date range to search, or un-check **Scheduled Procedure Step Start Date** and enter any other valid search term. Click **Next**. The **Add Patient - Step 3** guidance displays.

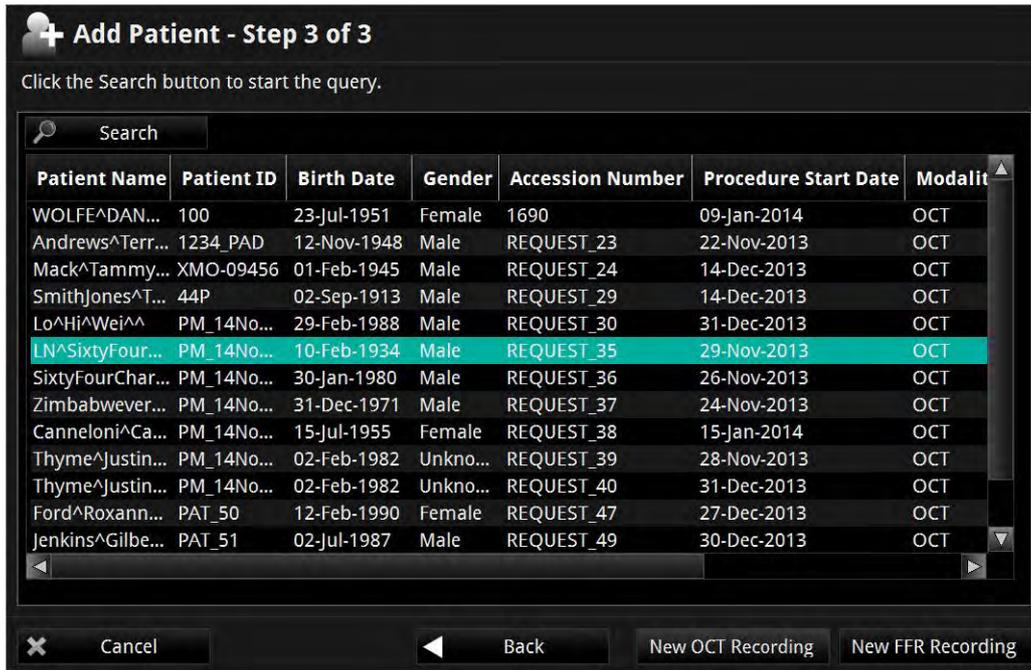
**+ Add Patient - Step 3 of 3**  
Click the Search button to start the query.

Search

Patient Name	Patient ID	Birth Date	Gender	Accession Number	Procedure Start Date	Modality
WOLFE^DAN...	100	23-Jul-1951	Female	1690	09-Jan-2014	OCT
Andrews^Terr...	1234_PAD	12-Nov-1948	Male	REQUEST_23	22-Nov-2013	OCT
Mack^Tammy...	XMO-09456	01-Feb-1945	Male	REQUEST_24	14-Dec-2013	OCT
SmithJones^T...	44P	02-Sep-1913	Male	REQUEST_29	14-Dec-2013	OCT
Lo^Hi^Wei^^	PM_14No...	29-Feb-1988	Male	REQUEST_30	31-Dec-2013	OCT
LN^SixtyFour...	PM_14No...	10-Feb-1934	Male	REQUEST_35	29-Nov-2013	OCT
SixtyFourChar...	PM_14No...	30-Jan-1980	Male	REQUEST_36	26-Nov-2013	OCT
Zimbabwever...	PM_14No...	31-Dec-1971	Male	REQUEST_37	24-Nov-2013	OCT
Canneloni^Ca...	PM_14No...	15-Jul-1955	Female	REQUEST_38	15-Jan-2014	OCT
Thyme^Justin...	PM_14No...	02-Feb-1982	Unkno...	REQUEST_39	28-Nov-2013	OCT
Thyme^Justin...	PM_14No...	02-Feb-1982	Unkno...	REQUEST_40	31-Dec-2013	OCT
Ford^Roxann...	PAT_50	12-Feb-1990	Female	REQUEST_47	27-Dec-2013	OCT
Jenkins^Gilbe...	PAT_51	02-Jul-1987	Male	REQUEST_49	30-Dec-2013	OCT

Cancel Back New OCT Recording New FFR Recording

4. Select the patient by clicking on the name. The **New OCT Recording** button and the **New FFR Recording** button become active.



5. Click either the **New OCT Recording** button or the **New FFR Recording** button to continue.

### From a DICOM Storage Server

To import patient information From a DICOM Storage Server, complete the following steps:

1. From the Select Patient menu, click the **Add Patient** button. The **Add Patient - Step 1** guidance displays.



Figure 8-11: Add Patient - Step 1 (Storage Server)

2. Select **DICOM Patient Info**. Click the **Next** button. The **Add Patient - Step 2** guidance displays.



3. If there is more than one DICOM Storage Server detected, select the desired server. Click the **Next** button. The **Add Patient - Step 3** guidance displays.



4. Enter a search term such as the **Patient Name** (shown), **Patient ID**, or leave blank to return the entire list of Patients. Click the **Search** button. The results display.



5. Click to highlight the patient, then click the **New OCT Recording** button or **New FFR Recording** button. The patient data is entered.



# Cleaning & Maintenance

---

9

**CAUTION:** Only a qualified service representative can service components of the system. Any attempt to open the system components by anyone other than a qualified St. Jude Medical service representative will void the warranty.

**CAUTION:** Do not perform cleaning or maintenance on the system in the patient environment.

## **Contacting St. Jude Medical Service**

Service can be contacted at:

E-mail: OCTservice@sjm.com  
Phone: +1 855 478 5833 US Toll-free  
+1 651 756 5833 International

## **Cleaning**

Cleaning of the ILUMIEN OPTIS System consists of:

- Cleaning system surfaces.
- Cleaning the DOC and its cable.

## **Routine Cleaning Procedure**

The OPTIS Integrated DOC, DOC Holster, and Table Side Controller should be cleaned following the facilities standard cleaning schedule, or at least every 30 days under normal use.

1. Turn off all system components with accessible power controls and unplug the power cables.
2. Clean system surfaces and the keyboard with a dry cloth, or a dry cloth slightly dampened with water.
3. Clean the monitor's LCD surface with a lint-free, non-abrasive cloth.

**CAUTION: DO NOT clean the LCD surface with detergents or other cleaning solutions.**

4. Clean exposed system cables with a soft cloth moistened with water or a mild detergent.
5. Clean the DOC and the DOC optical cable with a disinfectant wipe or Cidex<sup>TM</sup>1 (Glutaraldehyde 3.4%) and a soft cloth. Be particularly careful not to stress or sharply bend the DOC optical cable.

**NOTE:** Though enclosed in a bag during use, the DOC is the system component most exposed to dirt, fluids, and debris.

**NOTE:** The Table Side Controller is not water-proof. Be careful to not use excessive moisture when cleaning this device.

6. Clean all other exposed system cables with a soft cloth moistened with water or water and a mild detergent.

---

1. Cidex is a trademark of Johnson & Johnson Corporation.

## **Maintenance**

Maintenance of the system consists of:

- Cleaning the optical connection in the DOC and the Dragonfly™ Imaging Catheter.
- Replacing the optical adapter in the DOC.
- Inspecting exposed cable connections.
- Transferring log files.
- Identifying the installed software version.

## Optical Connection Cleaning Procedure

The optical connection between the DOC and the Dragonfly™ Imaging Catheter should be cleaned whenever there is a loss of image quality. Image quality should also be checked every 3 months, and the connection cleaned if there is a loss of quality.

**CAUTION:** This procedure should not be performed during a patient case on a sterile catheter.

**NOTE:** Before beginning this procedure ensure that you have the St. Jude Medical Optical Fiber Connector Cleaner available.

**NOTE:** Do not touch any of the optical connectors or the end of the Optical Fiber Connector Cleaner as this may damage them.

1. If a catheter is connected to the DOC, press the **Unload** button on the DOC and wait until the “lock” LED stops flashing. Once the “lock” LED is off, remove the catheter.
2. Click the **Menu** button and select the **Setup** option. The **Setup** dialog box displays.
3. Click the **Service** button. The **Service** menu displays.
4. In the **DOC Service** section of the **Service** menu, click the **Enter** button.

The DOC optical carriage moves all the way to the front and then locks into position for 1 minute. Rotation of the DOC rotary motor is stopped and the laser source is turned off.

5. Remove the sizing cap from the end of the Optical Fiber Connector Cleaner (see [Figure 9-1](#)).
6. Insert the cleaner into the center of the Optical Adapter in the DOC, ensuring it seats fully, and press until it clicks.

The cleaning material in the Optical Fiber Connector Cleaner is moved over the optical connection cleaning it.



Figure 9-1: Inserting Cleaner Into Optical Adapter

7. Remove the cleaner from the DOC.
8. In the DOC Service section of the **Service** menu, click the **Exit** button.

The DOC optical carriage moves all the way to the back and all DOC functions are returned to their normal state.

9. Click on **OK** or **Cancel** to close the **Setup** dialog box.
10. Open the sizing cap and place it on the end of the Optical Fiber Connector Cleaner (see [Figure 9-2](#)).
11. Insert the cleaner into the optical connection in the Dragonfly Imaging Catheter, ensuring it seats fully, and press until it clicks.

The cleaning material in the Optical Fiber Connector Cleaner is moved over the optical connection cleaning it.

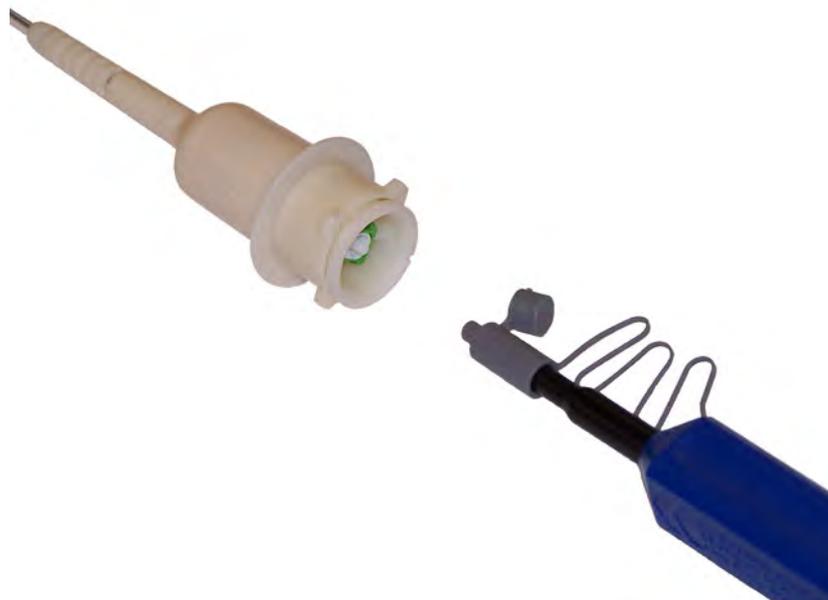


Figure 9-2: Inserting Cleaner Into Catheter

12. Remove the cleaner from the Dragonfly.
13. Reconnect the Dragonfly Imaging Catheter to the DOC.
14. Return the system to normal operation.

## **Optical Adapter Replacement Procedure**

The Optical Adapter should be replaced every 200 cycles or 1 year (which ever comes first), or if cleaning (see “[Optical Connection Cleaning Procedure](#)”) does not improve the image quality.

**CAUTION:** This procedure should not be performed during a patient case.

**NOTE:** Before beginning this procedure ensure that you have a replacement St. Jude Medical Optical Adapter and pliers supplied with it available.

**NOTE:** Do not touch any of the optical connectors as this may damage them.

1. If a catheter is connected to the DOC, press the **Unload** button on the DOC and wait until the “lock” LED stops flashing. Once the “lock” LED is off, remove the catheter.
2. Click the **Menu** button and select the **Setup** option. The **Setup** dialog box displays.
3. Click the **Service** button. The **Service** menu displays.
4. In the **DOC Service** section of the **Service** menu, click the **Enter** button.
5. Click on the **Enter DOC Service** button.

The DOC optical carriage moves all the way to the front and then locks into position for 1 minute. Rotation of the DOC rotary motor is stopped and the laser source is turned off.

6. Use pliers to grip either short surface of the adapter as shown (see [Figure 9-3](#)). Ensure the pliers seat over the lip on the end of the adapter.



Figure 9-3: Proper Gripping of Adapter for Removal

7. Pull the adapter straight out of the DOC.

**NOTE:** Ensure that the adapter does not rotate while it is being removed as this could damage the DOC or break the adapter.

8. Remove the cap and plug from the replacement Optical Adapter.

**NOTE:** Do not touch any of the optical connectors as this may damage them.

9. Using your fingers, align the new adapter with the DOC, ensuring the key is aligned with the slot in the optical carriage, and firmly press it into place.

There is a small click when the adapter seats.

**NOTE:** If the carriage moves when you press on the adapter, click the **Exit** button in the DOC Service section of the **Service** menu. Once it changes to Enter, click it again. Once the optical carriage is locked into position at the front of the DOC firmly press the adapter into place.

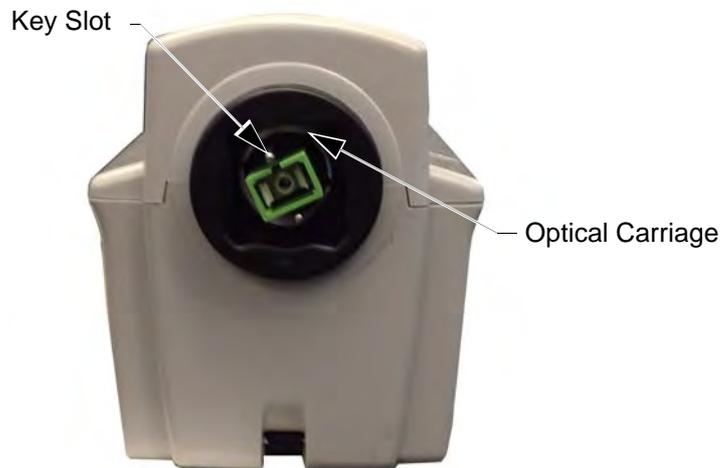


Figure 9-4: Alignment of Optical Adapter with Optical Carriage

10. In the DOC Service section of the **Service** menu, click the **Exit** button.

The DOC optical carriage moves all the way to the back and restores all DOC functions to their normal state.

11. Click on **OK** or **Cancel** to close the **Setup** dialog box.
12. Return the system to normal operation.

**Cable Connection Inspection Procedure**

1. Make sure the power cords are in good condition and are properly plugged in.
2. If a secondary ground connection is being used, ensure that both ends of the cable are secure.
3. Ensure all other connectors are fully seated and properly secured.

## Transferring Log Files

Log files may be transferred to a CD or external USB device so that you can archive them or E-mail them to St. Jude Medical Service.

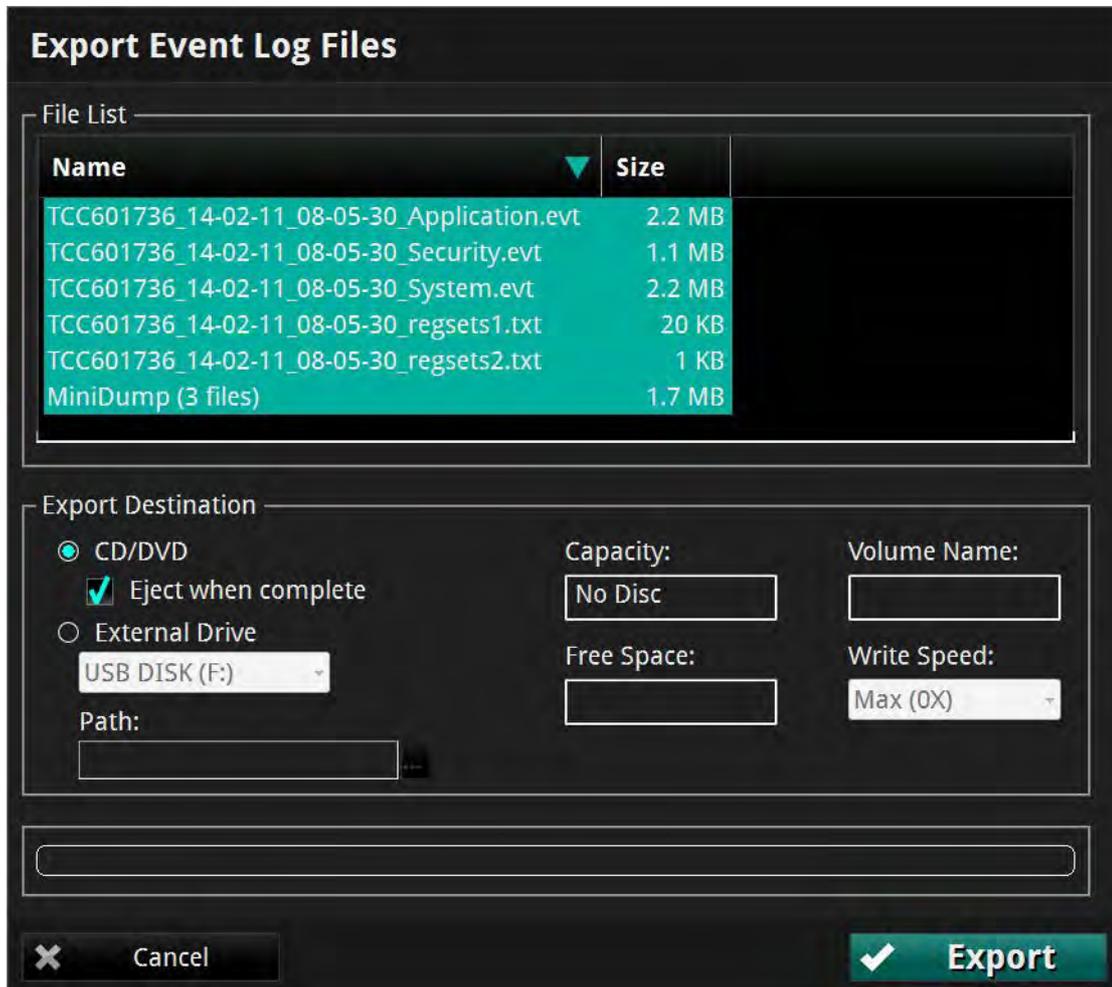


Figure 9-5: Transfer Event Log Files Menu

**NOTE:** When there are multiple repeated events of the same type the event logs are condensed into two messages. The first message logs the initial event, and the second message logs the event along with a repeat count.

1. Click the **Menu** button and select the **Setup** option. The **Setup** dialog box displays.
2. Click the **Service** button. The **Service** menu displays.
3. In the **Service Log** area, click the **Export** button.

The **Transfer Event Log Files** menu opens (Figure 9-5), which lists all the files to be exported.

4. Click to select a **Transfer Destination**:

**NOTE:** If an output device is not available, that option cannot be selected.

**NOTE:** Check that the **Free** space in the Drive Capacity section is large enough to save the files that are being exported.

- Click **CD/DVD** to export to a CD/DVD.

If the inserted CD/DVD is blank, you can edit the name in the **Volume Label** field.

- Click **External Drive** to export to an external USB device.

If necessary, click the **External Drive** selection box to select the correct drive location from the drop-down menu.

Click the **Ellipsis (...)** button to the right of the **Path** selection box to browse for a specific folder on the device.

5. Click the **Transfer** button to transfer the files or **Cancel** to cancel the operation.

## Identifying the Software Version

The **Startup** window (Figure 9-6) shows the version and copyright information for the system. If the system is running, type <Ctrl-S> to open the Startup window. The software version information might be required when contacting St. Jude Medical Service.



Figure 9-6: System Startup Window

## Infection Control

Follow the infection control procedures established in your institution for protection of both staff and patient.

Blood on system components, panels, and cables should be removed by using a gauze pad with soap and water, and drying with a soft cloth to prevent corrosion. The DOC cable can be cleaned with a disinfectant wipe or Cidex and a soft cloth.

## User Troubleshooting

Table 9-1 provides basic guidelines for troubleshooting the ILUMIEN OPTIS System. “Troubleshooting OCT Acquisition” on page 5-19 explains what to do if a catheter fails during pullback.

If your problem is not resolved after attempting the suggested remedies, contact your service representative or E-mail OCTservice@sjm.com.

Table 9-1: User Troubleshooting Tips

Symptom	Possible Causes	Remedy
<b>General</b>		
Screen blank, power indicator on monitor not lit.	Display not turned on.	Press the power button on the monitor to turn on monitor power.
	Display power cord unplugged.	Plug the monitor power cord into the back of monitor.
	System power not turned on.	Turn on the system power with the main power switch, located next to the power cord connection.
	System power cord not plugged in or not tight at either system end or wall connection.	Check to make sure plug is tightly connected to both the system and to the wall outlet.
	Main system power fuse is blown.	Refer to the Electrical label on the system (located next to the power cord connection) for fuse information. Contact your service representative for instructions.
Screen blank, power indicator on monitor lit.	Outlet power disrupted.	Check voltage at the wall outlet.
	PC auto-boot failed.	Turn off the main power switch and wait fifteen seconds. Turn the main power switch back on, and press the On/Standby button on the right side of the keyboard to start the system. If the system still does not start, contact your St. Jude Medical service representative for instructions.
	Monitor not enabled.	Click <Alt+C> on the keyboard to enable the display. If the screen remains blank, turn off the main power switch and wait fifteen seconds. Turn the main power switch back on, and press the On/Standby button on the right side of the keyboard to start the system.

Table 9-1: User Troubleshooting Tips (*continued*)

<b>Symptom</b>	<b>Possible Causes</b>	<b>Remedy</b>
<b>Connections</b>		
During data export, the system does not list the USB drive plugged into the USB port.	The connected USB drive is not compatible with the system, or the USB drive requires formatting.	Connect only a USB Drive supported by Windows 7 SP1. Refer to the instructions for use that came with your USB drive to determine if formatting is required.
Screen message “Imaging engine initialization failed” is displayed at startup.	This message can be caused by several problems, including loose or damaged system connections.	Shut down the system, turn off the main system power, and wait 15 seconds. Then turn the system back on. If the error is displayed again, note the error code displayed (if any) then contact your St. Jude Medical service representative for instructions.
<b>DOC</b>		
Optical fiber does not rotate when <b>Live View</b> is pressed.	The <b>Stop</b> button on the DOC was pressed.  Imaging catheter defective, optical fiber does not rotate.	Check screen for message and follow instructions.  Replace imaging catheter.
DOC makes excessive noise without imaging catheter connected.	DOC mechanism failure.	Contact your service representative to obtain a replacement DOC.
<b>Imaging</b>		
OCT image dim, with no background noise visible.	Monitor contrast and brightness set incorrectly.  Image contrast levels set incorrectly.	Set monitor contrast and brightness using monitor controls on the monitor.  Check the <b>Presentation Settings</b> on the <b>Settings</b> tab of the <b>OCT Settings</b> menu (normal settings are: Black level = 5%, White level = 90%. With the optical fiber rotating (scanning mode), reduce the Black level until background noise just becomes visible. Lack of background noise during optical fiber rotation indicates a defective imaging engine. See other possible causes and remedies below.
	Defective imaging catheter causing system saturation.	Remove imaging catheter from DOC. If background noise appears, the imaging catheter is defective. Replace catheter.

## Cleaning & Maintenance User Troubleshooting

---

Table 9-1: User Troubleshooting Tips (*continued*)

Symptom	Possible Causes	Remedy
	Dirty connection between DOC and imaging catheter.	Refer to <a href="#">“Optical Connection Cleaning Procedure”</a> on page 9-5 to clean the connection.

## System Disposal



Disposal of the equipment must be in accordance with local laws.



## Setup Dialog Box and Submenus

The **Setup** dialog box provides access to all of the configuration and service functions for the system. The **Setup** dialog box is accessed by selecting the **Setup** option from the **Menu** button options on the **Select Patient** or **Patient Summary** screen.

**NOTE:** When creating or reviewing an OCT recording, the **Settings** button opens a context-sensitive OCT **Settings** menu. When creating an FFR recording, the **Settings** button opens the **FFR Settings** menu.

All menus within the **Setup** dialog box share the following common buttons located at the bottom of the dialog box.

Table 10-1: Setup Dialog Box Common Options

<b>Refresh</b>	Click to update the values displayed in the current tab.
<b>OK</b>	Click to approve any changes and close the <b>Setup</b> dialog box.
<b>Cancel</b>	Click to cancel any changes and close the <b>Setup</b> dialog box.
<b>Apply</b>	Click to approve any changes and keep the <b>Setup</b> dialog box open.

## Setup - Acquisition Menu

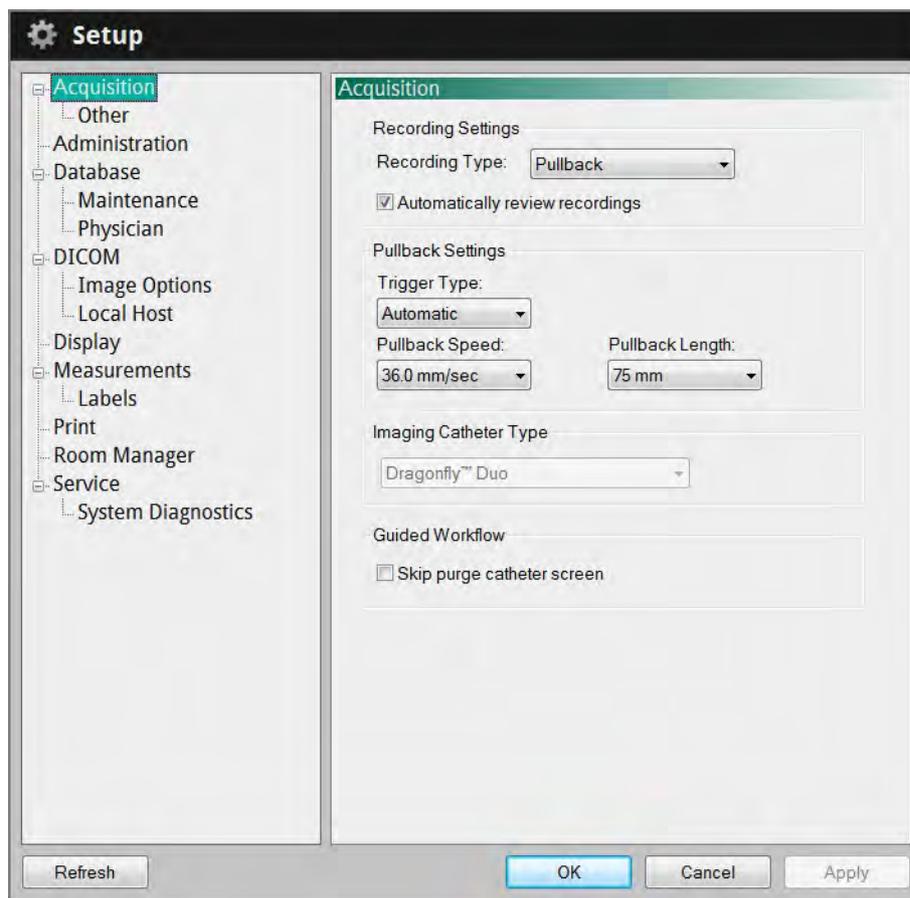


Figure 10-1: Setup - Acquisition Menu

Table 10-2: Setup - Acquisition Menu Settings

<p><b>Recording Type</b></p>	<p>Sets the recording type for image acquisition:</p> <p><b>Pullback</b> - Performs a recording while the imaging core of the catheter is pulled back within the catheter sheath.</p> <p><b>Stationary</b> - The system records the live view image for 6 seconds without pulling the imaging core of the catheter back.</p>
<p><b>Automatically review recordings</b></p>	<p>Turns on and off automatic review after recording real-time images.</p>

Table 10-2: Setup - Acquisition Menu Settings (*continued*)

<p><b>Trigger Type</b></p>	<p>Sets the Pullback trigger type to use for image acquisition:</p> <p><b>Automatic</b> - The system automatically begins a recording when it detects that the vessel has been cleared by the flush injection.</p> <p><b>Manual</b> - The System does not perform a pullback until you click the <b>Start Pullback</b> button on the screen, or press the <b>Enable</b> button on the DOC.</p>
<p><b>Pullback Speed</b></p>	<p>Sets the pullback speed.</p> <p><b>NOTE:</b> When connected to an original C7 Dragonfly imaging catheter, the <b>Pullback Speed</b> settings are limited to <b>10.0 mm/sec, 20.0 mm/sec, and 25.0 mm/sec</b>. For Dragonfly Duo or OPTIS, available speeds are <b>18</b> and <b>36 mm/sec</b>.</p>
<p><b>Pullback Length</b></p>	<p>Sets the length of the pullback.</p> <p><b>NOTE:</b> When connected to an original C7 Dragonfly imaging catheter, the <b>Pullback Length</b> setting is limited to <b>54 mm</b>. For Dragonfly Duo or OPTIS <b>Pullback Length</b> is <b>54 mm</b> when speed is <b>18 mm/sec</b> and <b>75 mm</b> when speed is <b>36 mm/sec</b>.</p>
<p><b>Imaging Catheter Type</b></p>	<p><b>NOTE:</b> This displays the type of catheter when a catheter has been connected.</p> <p><b>NOTE:</b> Only Dragonfly imaging catheters are approved for use with the ILUMIEN OPTIS System for cardiovascular imaging.</p>
<p><b>Guided Workflow</b></p>	<p>Turns on and off the appearance of the purge catheter guidance message.</p>

## Setup - Acquisition/Other Menu

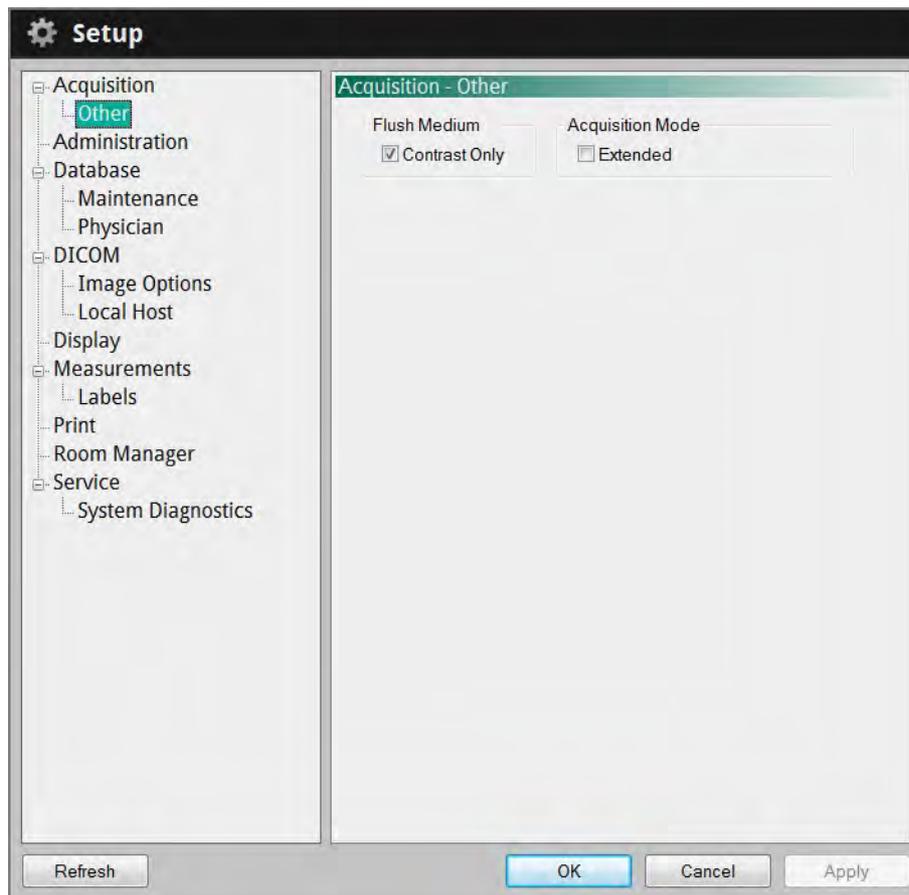


Figure 10-2: Setup - Acquisition/Other Menu

Table 10-3: Setup - Acquisition/Other Menu Settings

<p><b>Flush Medium</b></p>	<p>Controls the contrast options available in the <b>Flush Medium</b> drop-down box.</p> <p><b>Checked</b> - Only 100% contrast is available, and the <b>Flush Medium</b> drop-down box under <b>Settings</b> is unavailable.</p> <p><b>Unchecked</b> - All configured contrast options are available for selection in the <b>Flush Medium</b> drop-down box.</p> <p><b>NOTE:</b> To change the <b>Flush Medium</b> setting before performing an OCT recording, see <a href="#">“Confirm Recording Settings”</a> on page 5-12.</p>
<p><b>Acquisition Mode</b></p>	<p>This option should remain unchecked unless instructed otherwise by SJM personnel.</p>

## Setup - Administration Menu

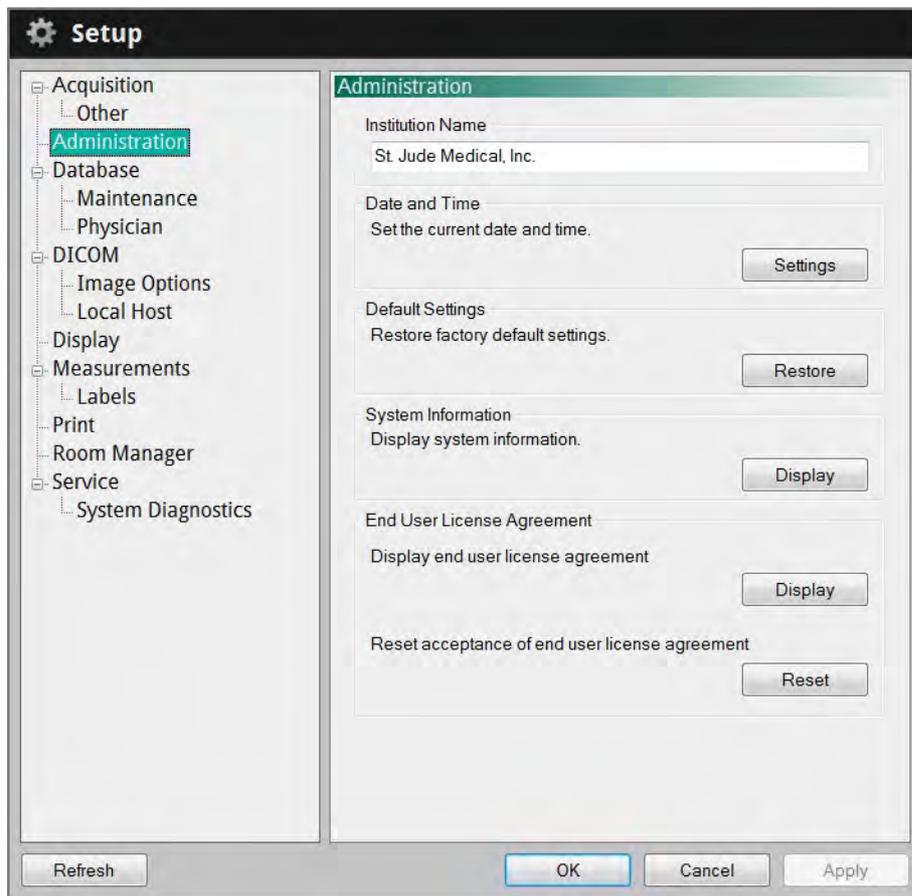


Figure 10-3: Setup - Administration Menu

Table 10-4: Setup - Administration Menu Settings

<b>Institution Name</b>	Click here and enter or clear the institution name.
<b>Date and Time</b>	Opens the system <b>Date and Time Properties</b> menu.
<b>Default Settings</b>	Resets all user-entered configuration values except the date and time to the original factory default values. <b>CAUTION: Restoring factory default settings resets ALL user-entered configuration values except the date and time. This button should be used only by qualified service personnel or under their direction.</b>

Table 10-4: Setup - Administration Menu Settings (*continued*)

<b>System Information</b>	Displays current date, system name, institution name, and software version
<b>End User License Agreement</b>	Displays EULA. Resets EULA. The EULA must be accepted upon first operation after reset in order to access the system.

## Setup - Database Menu

The **Database** menu in the **Setup** dialog box shows a list of all files in the database. From this menu, you can select patients or individual recordings for export, delete recordings, or select a patient for editing.

Click on individual recordings to select them, or click on a patient's name to select all recordings for that patient, or click on the database to select all recordings.

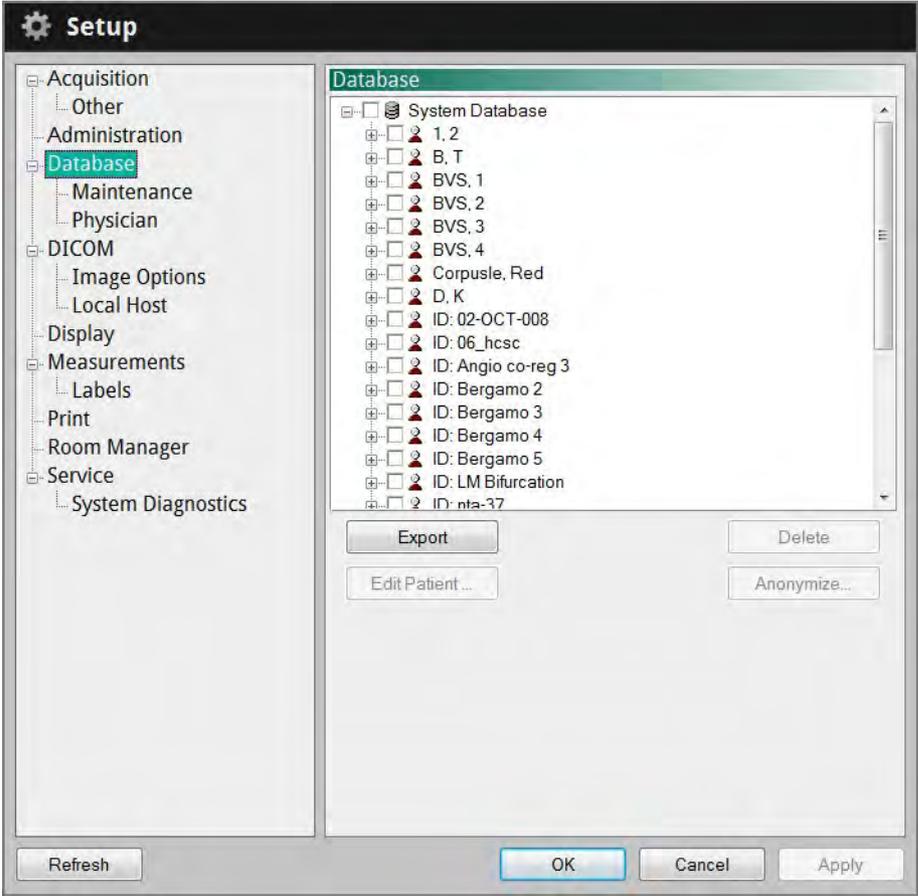


Figure 10-4: Setup - Database Menu

Table 10-5: Setup - Database Menu Settings

<b>Database display</b>	The database display shows all recordings, grouped by patient name, and then by the date of the recordings.
<b>Export</b>	Exports the selected recordings.
<b>Delete</b>	Deletes the selected recordings.

Table 10-5: Setup - Database Menu Settings (*continued*)

<b>Edit Patient...</b>	Click on a patient's name to edit the patient's information. See <a href="#">“Editing Patient Information”</a> on page 3-8 for more information.
<b>Anonymize...</b>	Click on a patient's name to remove the patient's identifying information from the record. See <a href="#">“Exporting Files in Native (Raw) Format”, Step 3</a> , on page 8-8, for an explanation of the anonymization function.

## Setup - Database/Maintenance Menu

You can check the status of the system’s OCT database using the **Database Maintenance** menu in the **Setup** dialog box. From this menu, you can check the size of the database file, as well as the number of files referenced by the database and the combined size of these files. This menu also warns you of possible issues found with files in the system’s database.

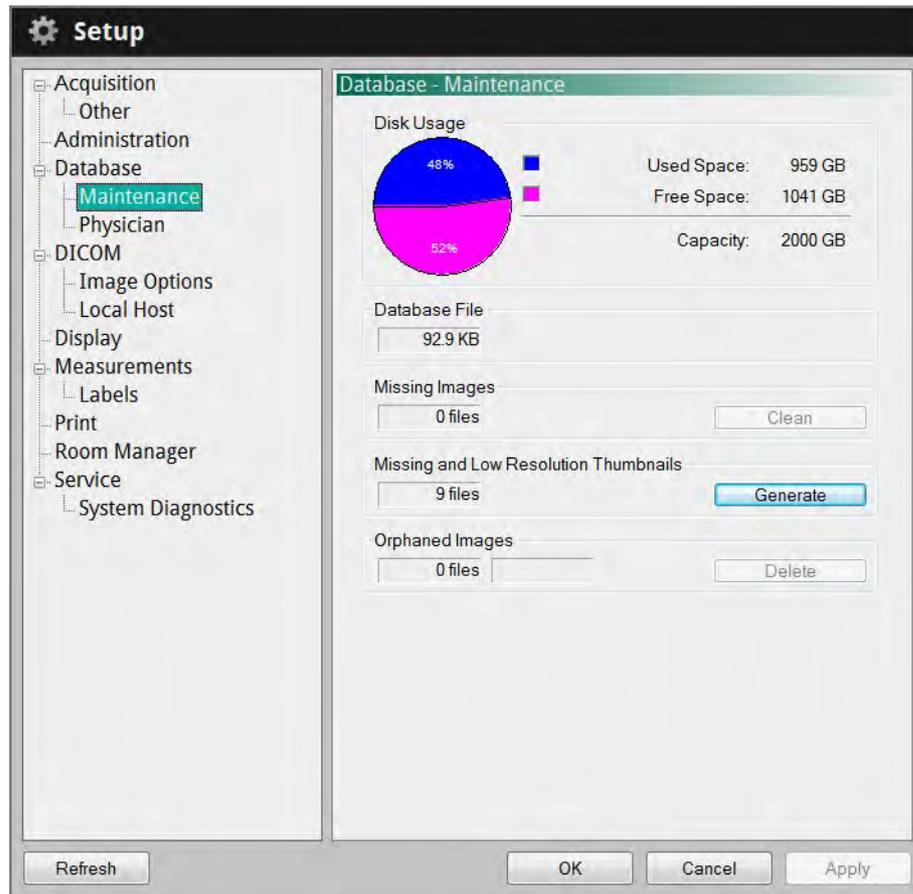


Figure 10-5: Setup - Database/Maintenance Menu

Table 10-6: Setup - Database/Maintenance Menu Settings

<b>Disk Usage</b>	Provides both a graphic and text showing the amount of free and used space on the system’s hard disk.
-------------------	---

**User Interface Reference**  
**Setup - Database/Maintenance Menu**

---

Table 10-6: Setup - Database/Maintenance Menu Settings (*continued*)

<b>Database File</b>	Indicates the size of the database file.
<b>Missing Images</b>	OCT image files that are referenced by the database but could not be found. Click the <b>Clean</b> button to remove references to files that cannot be found.
<b>Missing and Low Resolution Thumbnails</b>	Thumbnail image files that are missing or are low resolution. Click the <b>Generate</b> button to generate new thumbnail image files.
<b>Orphaned Images</b>	OCT image files that were found on the system but are not referenced by the database. Click the <b>Delete</b> button to delete these orphaned files to free up additional hard disk space.

## Setup - Database/Physician Menu

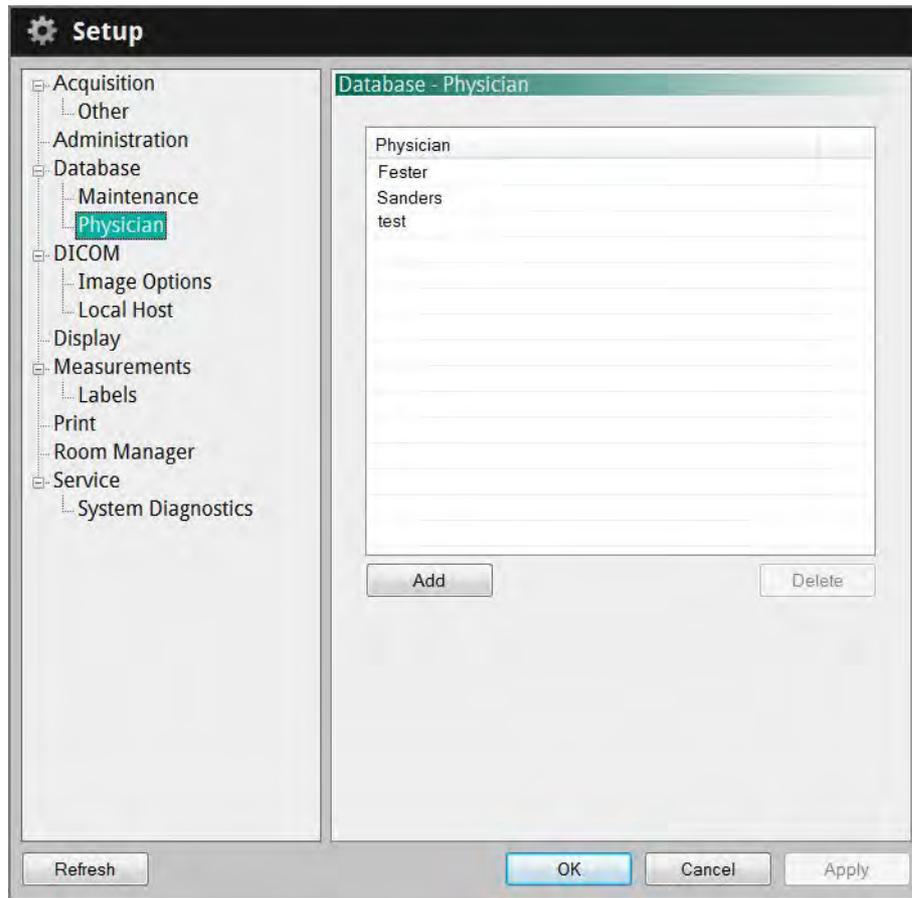


Figure 10-6: Setup - Database/Physician Menu

Table 10-7: Setup - Database/Physician Settings

<b>Physician list</b>	Shows a list of all physician names entered in the system. <ul style="list-style-type: none"><li>• Click <b>Add</b> to add a new physician name.</li><li>• Click a name and click <b>Delete</b> to remove the name.</li></ul>
-----------------------	---

## Setup - DICOM Menu

The **DICOM** menus in the **Setup** dialog box are used to configure the network settings of the system and the Remote DICOM Server so that the system can initiate a network connection between the two systems. If multiple network adapters exist on the system (not a St. Jude Medical standard configuration) only the primary network adapter can be used.

**CAUTION:** The **ILUMIEN OPTIS System** should only be connected to a secure intranet. Direct connection to the Internet may interfere with correct operation and/or result in inappropriate access to patient information, and voids system warranties.

Please note St. Jude Medical makes no representation or warranty that use of the **ILUMIEN OPTIS System** complies with applicable privacy, security and confidentiality laws, but encourages you to assess your own risk as you use, disclose, control, process or transfer patient health information with the system.

It is strongly recommended you contact your IT department to set up the system on the DICOM server. Refer to the *St. Jude Medical DICOM Conformance Statement* for more information on configuration.

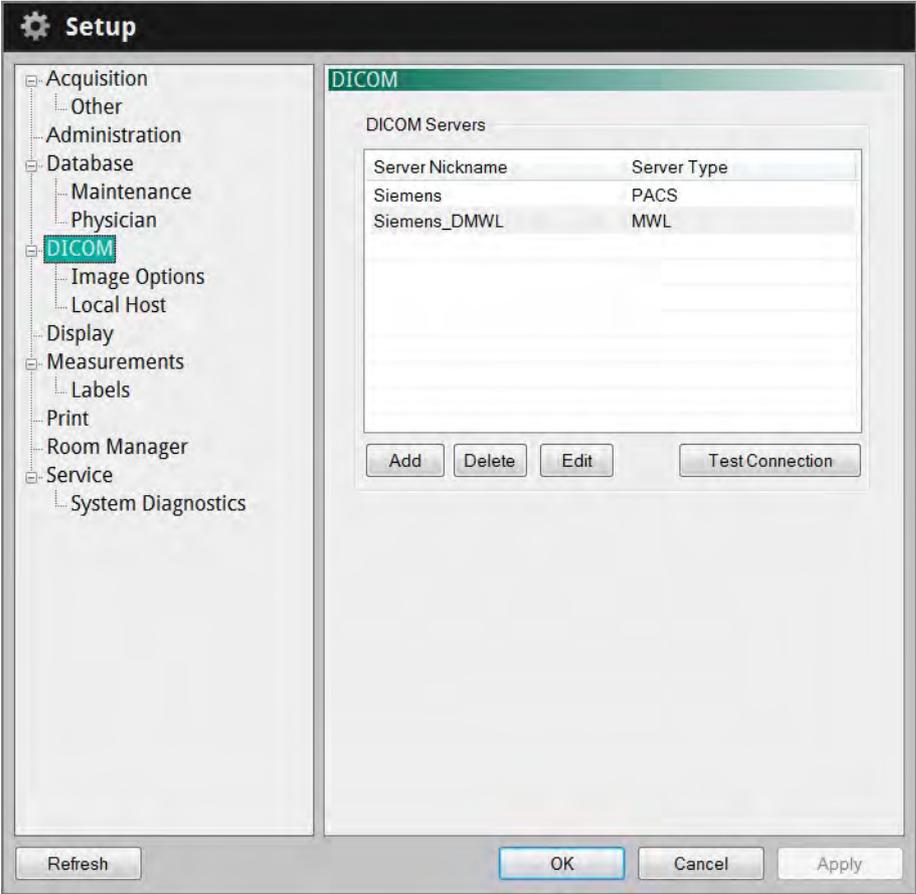


Figure 10-7: Setup - DICOM Menu

Table 10-8: Setup - DICOM Menu Settings

<b>Server Nickname</b>	Displays the Nickname for the server.
<b>Server Type</b>	Displays the type of DICOM server; either <b>PACS</b> or Modality Worklist ( <b>MWL</b> )
<b>Add</b>	Click to <b>Add</b> a server.
<b>Delete</b>	Click to <b>Delete</b> a server.
<b>Edit</b>	Click <b>Edit</b> to display the <b>Configure DICOM</b> menu.
<b>Test Connection</b>	Click to send a <b>Test</b> signal to the server and confirm that you have successful communication.

**Configure DICOM**

Remote Host

Hostname

Network Timeout(secs):

15

Specify IP Address

Ping Host

DICOM Store Settings

Local AE Title

Remote AE Title

Remote Port

104

Remote Nickname

Response Timeout(secs):

600

Query Server

OK Cancel

Figure 10-8: Configure DICOM Menu

Table 10-9: Setup - Configure DICOM Menu Settings

<b>Hostname</b>	The host name of the network server that contains the Remote DICOM Server. Use this option to identify the server if the network supports the Domain Name System (DNS) and the host name is known; otherwise, use <b>Specify IP Address</b> (see below).
<b>Network Timeout (secs)</b>	The maximum time to allow for a network ping response from the network server that contains the Remote DICOM Server. The minimum setting is 1 second while the maximum is 120 seconds. The default is 15 seconds.
<b>Specify IP Address</b>	The IP address of the network server that contains the Remote DICOM server. If the Remote Hostname is used to identify the server, this field will be automatically filled in if the Ping Host button (Figure 10-8) is successfully used to verify the network connection. Check the Specify checkbox to the right of this option to explicitly specify the remote IP address.
<b>Ping Host</b>	This button can be used to test the network connection between the System and the Remote DICOM Server Host. If Hostname or IP address is not correct, a <b>Ping Failed</b> error message displays.
<b>Local AE Title</b>	The AE (Application Entity) title used to identify the Local DICOM Store SCU (Service Class User) used by the ILUMIEN OPTIS System.
<b>Remote AE Title</b>	The AE title used to identify the Remote DICOM Server SCP (Service Class Provider) to which the system connects.
<b>Remote Port</b>	The port number on which the Remote DICOM Server Host will be listening for connection requests. The default is 104.
<b>Remote Nickname</b>	The nickname of the remove server defined by the user (16 characters maximum).
<b>Query Server</b>	Click to send a request from the system to the remote server to query the presentation context supported by the remote server. A test connection request is also sent from the server along with the query.

Table 10-9: Setup - Configure DICOM Menu Settings (*continued*)

<b>Response Timeout (secs)</b>	<p>The maximum time to allow for a response from the Remote DICOM Server after sending a DICOM request. The minimum setting is 15 seconds while the maximum is 1800. The default is 600 seconds.</p> <p>Press the <b>Test Connection</b> button to test the connection between the Local DICOM Store SCU and the Remote DICOM Server SCP. For a successful test:</p> <ul style="list-style-type: none"><li>• A TCP/IP connection is successfully made from the system to the remote server using the server IP address and port.</li><li>• A DICOM Associate Connection between the Local SCU and the Remote SCP is accepted by the SCP.</li><li>• A C-Echo request (Verification class) from the Local SCU results in a successful response from the Remote SCP.</li></ul> <p><b>NOTE:</b> If the DICOM server does not support Multi-Frame True Color Secondary Capture, a warning message appears, telling you to use SC-Image-Storage SOP class instead. If the DICOM server supports neither the Multi-Frame True Color Secondary Capture nor Secondary Capture, a warning message appears: "The presentation syntax supported by the remote server is not compatible. Images cannot be exported to the remote DICOM server."</p> <p><b>NOTE:</b> The <b>Test Connection</b> button is unavailable if the <b>Obtain an IP Address Automatically</b> checkbox is checked. See "<a href="#">Setup - DICOM/Local Host Menu</a>" on page 10-19.</p>
--------------------------------	--

Setup - DICOM/Image Options Menu

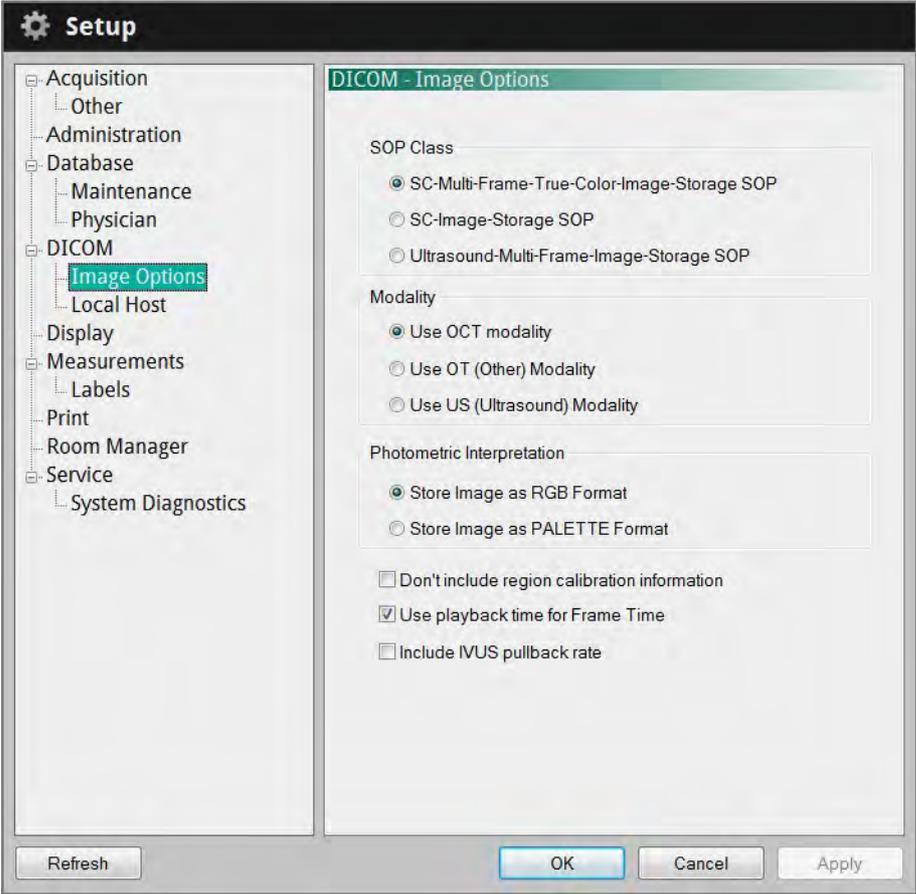


Figure 10-9: Setup - DICOM/Image Options Menu

Table 10-10: Setup - DICOM/Image Options Menu Settings

<p><b>SOP Class</b></p>	<p><b>SC-Multi-Frame-True- Color-Image-Storage SOP:</b> When this button is selected, DICOM images are exported as Multi-Frame True Color Secondary Capture. This selection is checked by default.</p> <p><b>SC-Image-Storage SOP:</b> When this button is selected, DICOM images are exported as Single-Frame Secondary Capture.</p> <p><b>Ultrasound-Multi-Frame-Image-Storage SOP:</b> When this button is selected, the images are exported as Ultrasound Multi-Frame DICOM data sets.</p>
-------------------------	--

**User Interface Reference**  
**Setup - DICOM/Image Options Menu**

---

Table 10-10: Setup - DICOM/Image Options Menu Settings (*continued*)

<b>Modality</b>	Check to select the image modality: OCT, OT (Other), or US (Ultrasound). The default DICOM export modality is OCT.
<b>Photometric Interpretation</b>	Check to select the export format for DICOM images: RGB or PALETTE. The default format is RGB. PALETTE yields smaller file sizes.
<b>Don't include region calibration information</b>	If this box is checked, region calibration information is not included in the DICOM data.

## Setup - DICOM/Local Host Menu

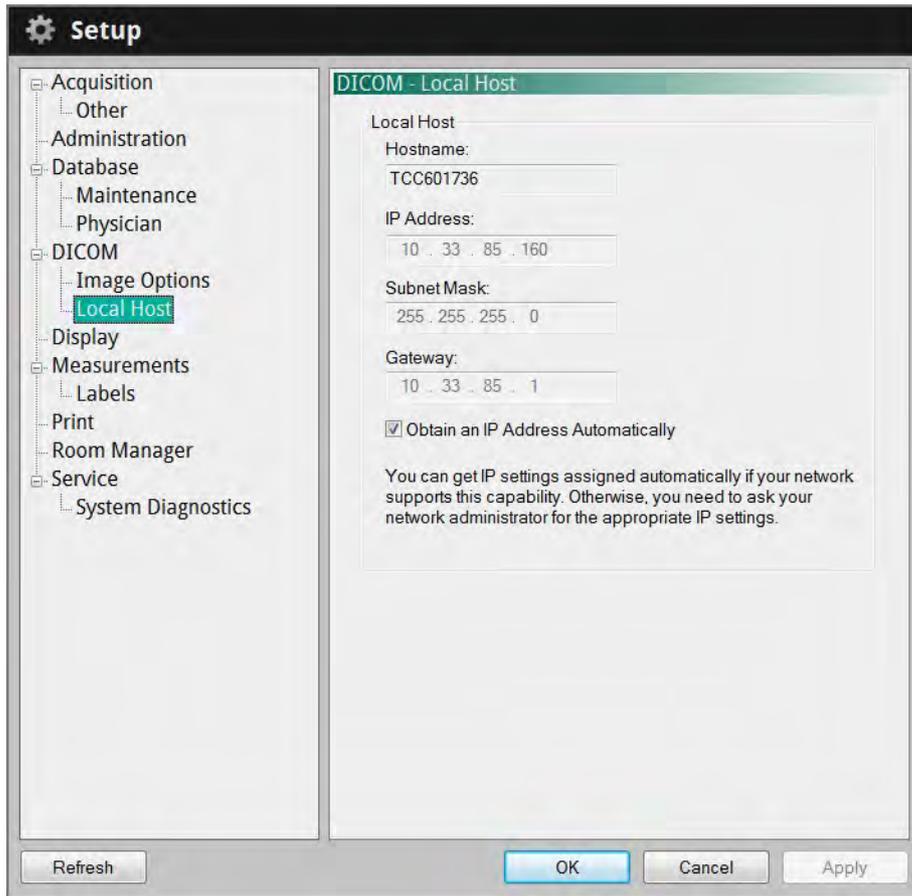


Figure 10-10: Setup - DICOM/Local Host Menu

Table 10-11: Setup - DICOM/Local Host Menu Settings

<b>Hostname</b>	The name used to identify the ILUMIEN OPTIS System on the network. This name is shown for display purposes only and cannot be modified.
<b>IP Address</b>	The IP address of this computer on the network. The default value is determined using DHCP when the network adapter is initialized.  <b>NOTE:</b> If the <b>Obtain an IP Address Automatically</b> option is checked, this is unavailable.  <b>NOTE:</b> Setting the IP address to an address used by another system on the network may cause network instability.

Table 10-11: Setup - DICOM/Local Host Menu Settings (*continued*)

<b>Subnet Mask</b>	The subnet mask number which is combined with the Local IP Address to identify which network segment this computer is on. The default value is determined using DHCP when the network adapter is initialized. <b>NOTE:</b> If the <b>Obtain an IP Address Automatically</b> option is checked, this is unavailable.
<b>Gateway</b>	The IP address of the default gateway, IP router, to be used to forward network traffic beyond the local network. The default value is determined using DHCP when the network adapter is initialized. <b>NOTE:</b> If the <b>Obtain an IP Address Automatically</b> option is checked, this is unavailable.
<b>Obtain an IP Address Automatically</b>	When checked, the <b>IP Address</b> of the Local Host, <b>Subnet Mask</b> , and <b>Gateway</b> are obtained by the system using DHCP, and cannot be edited on the ILUMIEN OPTIS System. Uncheck this option to modify the <b>IP Address</b> of the Local Host, <b>Subnet Mask</b> , and <b>Gateway</b> .

Setup - Display Menu

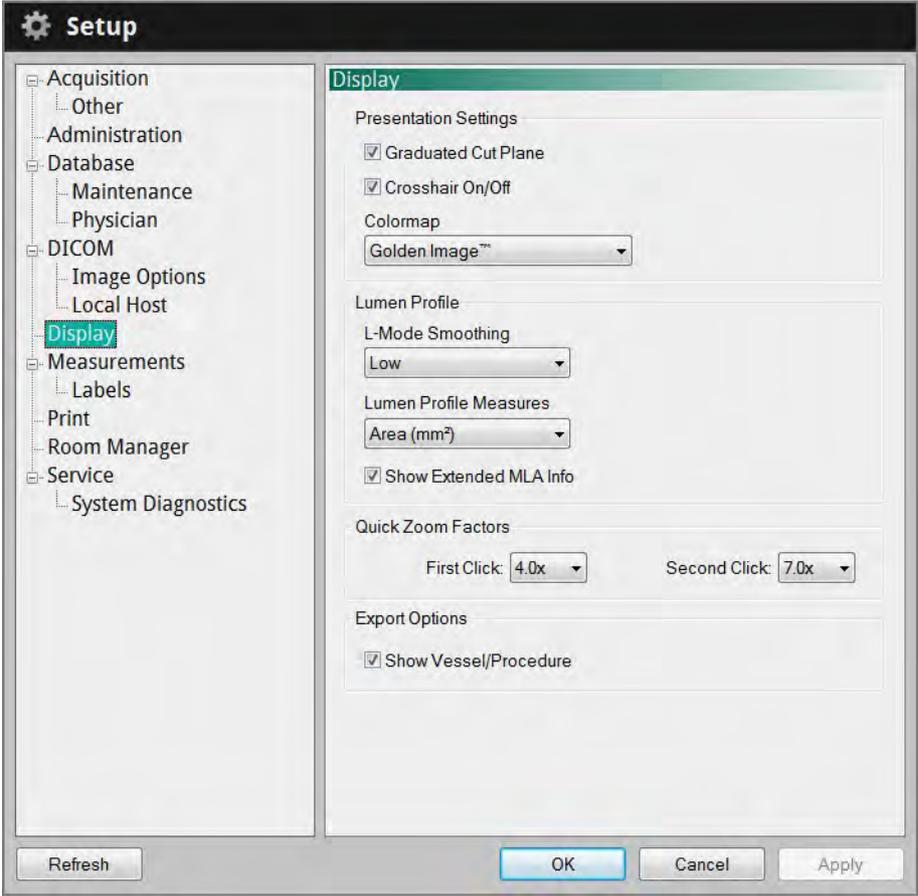


Figure 10-11: Setup - Display Menu

Table 10-12: Setup - Display Menu Settings

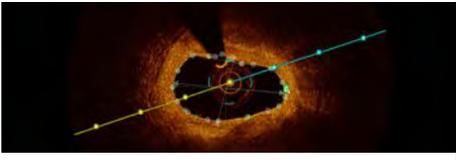
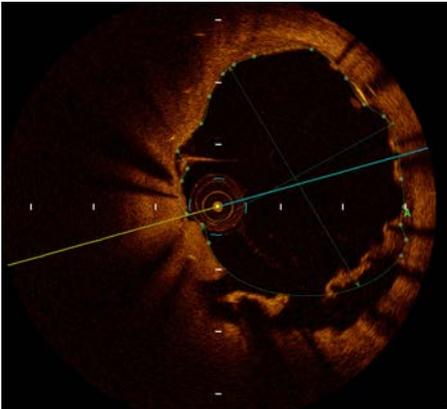
<p><b>Graduated Cut Plane</b></p>	<p>Check this checkbox to add graduated marks to the cut plane indicator to use as scale marks.</p> 
-----------------------------------	--

Table 10-12: Setup - Display Menu Settings (*continued*)

<p><b>Crosshair On/Off</b></p>	<p>Check this checkbox to add scaled crosshairs to the cross section view.</p> 
<p><b>Colormap</b></p>	<p>Click the arrow on the <b>Colormap</b> drop-down menu to display the list of colors. Click a color to select it.</p>
<p><b>L-Mode (Smoothing)</b></p>	<p>Click to select the amount of smoothing (averaging) for L-Mode views. The default setting is <b>Low</b>.</p>
<p><b>Lumen Profile Measures</b></p>	<p>Select the type of measurements to be represented in the Lumen Profile view. The default value is Mean Diameter. When Mean Diameter is selected, %DS values are displayed in the Lumen Profile view. These values are changed to %AS when Area is selected.</p>
<p><b>Show Extended MLA Info</b></p>	<p>Check to include the %DS (or %AS) values of the MLA with respect to both the Proximal and Distal Reference frames. This is in addition to the %DS (or %AS) value of the MLA with respect to the mean of the Proximal and Distal reference frames.</p>
<p><b>Quick Zoom Factors</b></p>	<p>Select the first click and second click zoom factor for the 2D, 3D, and Angio displays. Available choices for the first click setting are: 2.0x; 3.0x; 4.0x; 5.0x; 6.0x; and 7.0x. The second click setting can be set to a zoom factor greater than or equal to the first click setting. If the first click and second click settings are set to the same value then the first click will zoom in and the second click will zoom out. If the two settings are different then the first and second clicks will zoom in to their zoom factor, and the third click will zoom out. The default values are 3.0x and 6.0x respectively.</p>
<p><b>Show Vessel/Procedure</b></p>	<p>Check this box to present the Vessel and Procedure value on the exported DICOM image. The default is unchecked.</p>

## Setup - Measurements Menu

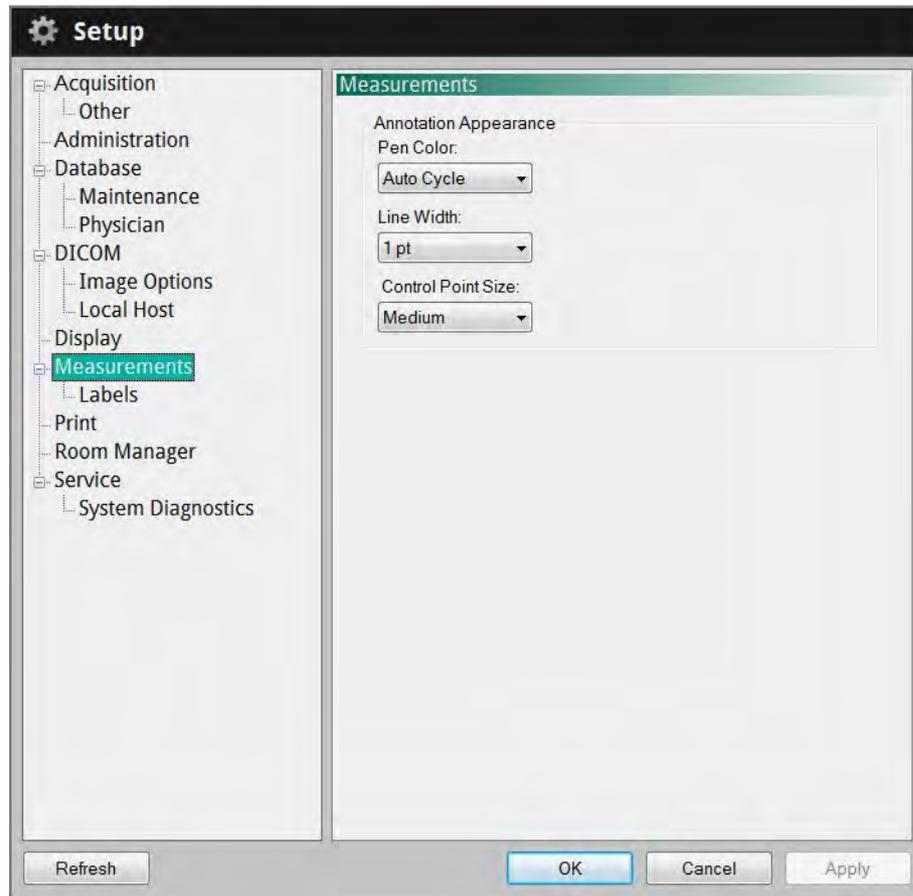


Figure 10-12: Setup - Measurements Menu

Table 10-13: Setup - Measurements Menu Settings

<b>Pen Color</b>	List of pen colors that can be used for subsequent measurements and calculations. If <b>Auto Cycle</b> is selected, the pen color used for measurements is automatically changed to the next color after a new measurement is completed.
<b>Line Width</b>	List of line widths that can be used for subsequent measurements and calculations. The default is <b>1 pt</b> .
<b>Control Point Size</b>	List of point sizes that can be used for subsequent length and two-point circular area measurements. The default is <b>Medium</b> .

## Setup - Measurements/Labels Menu

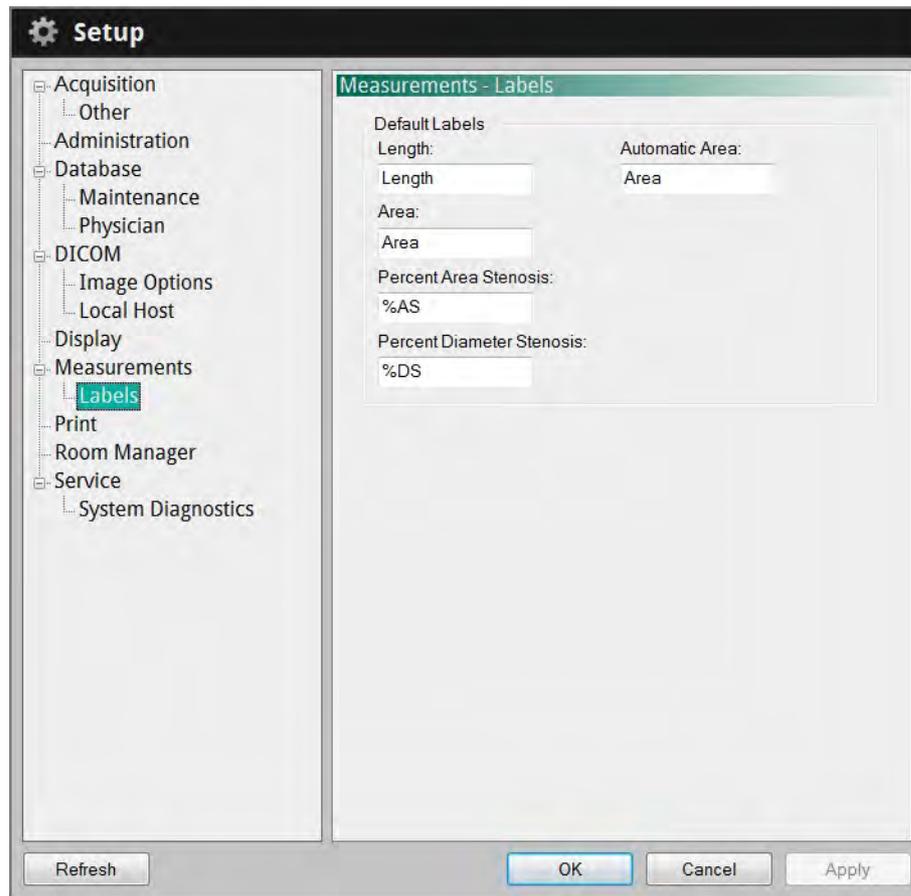


Figure 10-13: Setup - Measurements/Labels Menu

Table 10-14: Setup - Measurements/Labels Menu Settings

<b>Default Labels</b>	Label text that will be used when making the indicated type of measurement. Click and type in any field to change that label.
-----------------------	--

### Setup - Print Menu

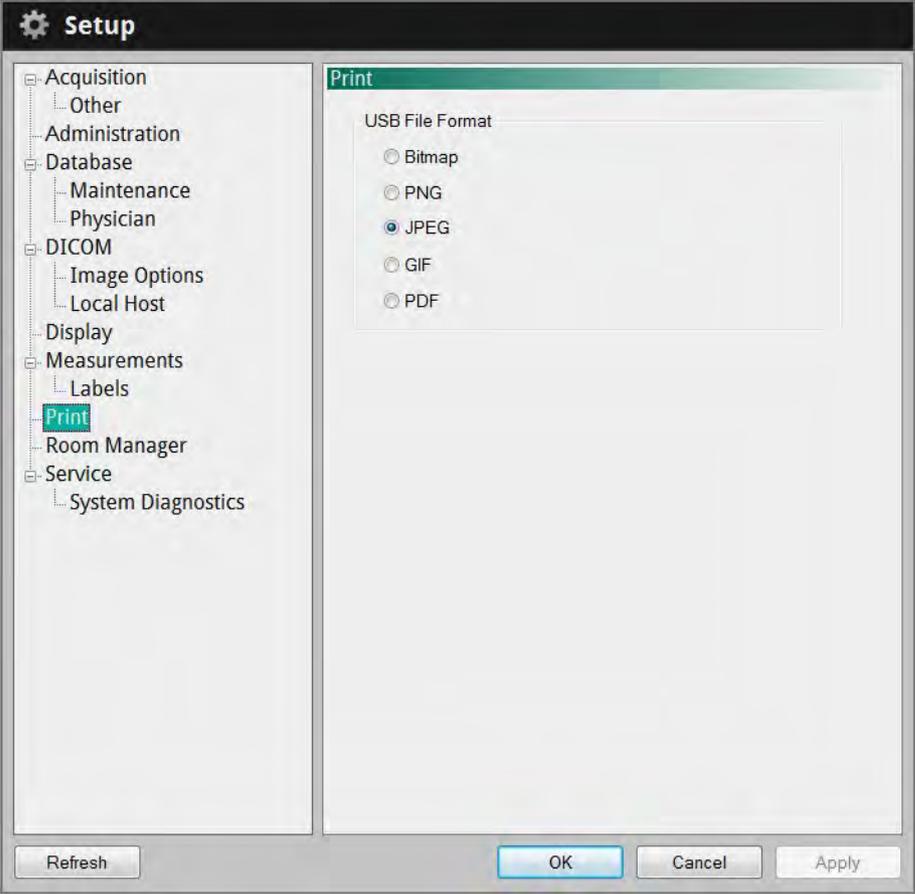


Figure 10-14: Setup - Print Menu

Table 10-15: Setup - Print Menu Settings

<b>USB File Format</b>	Click to choose the format used when printing a file to a USB drive.
------------------------	--

## Setup - Service Menu

The **Service** menu in the **Setup** dialog box displays service logs to help St. Jude Medical Service analyze problems which may occur during system operation. If you encounter problems while using the system, you may be asked to access and report information from this menu.

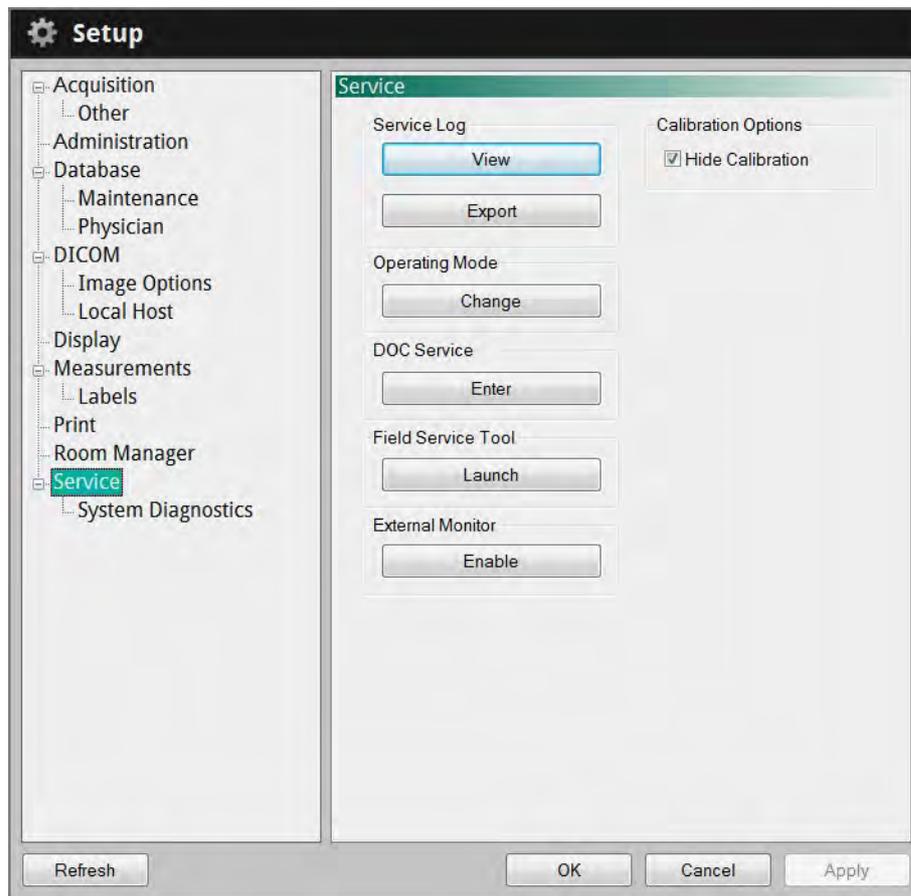


Figure 10-15: Setup - Service Menu

Table 10-16: Setup - Service Menu Settings


<p><b>Service Log (View)</b></p>	<p>Click the <b>View</b> button to open the Service Event Viewer.</p> <p><b>Buttons</b> →</p> <p><b>Headers</b> →</p>  <p>Click on a column header to sort the patients according to data in that column. Click the buttons at the top to show or hide categories of information:</p> <ul style="list-style-type: none"> <li>• <b>Errors</b> - lists all the operating errors logged by the system.</li> <li>• <b>Warnings</b> - lists all the Warning prompts displayed during operation.</li> <li>• <b>Information</b> - lists selected system status events which may help diagnose problems.</li> <li>• <b>No Category</b> - Displays entries that do not have a Category listed.</li> <li>• <b>Patient</b> - lists all changes and import/export of patient data.</li> </ul>
<p><b>Service Log (Export)</b></p>	<p>Click the <b>Export</b> button to open the <b>Transfer Event Log Files</b> menu, to export the Event Log files to a specified drive. See <a href="#">“Transfer Event Log Files Menu” on page 9-11</a> for more information.</p>
<p><b>Operating Mode</b></p>	<p>Click the <b>Change</b> button to open the <b>System Configuration Utility</b> window to change the application type and/or operating mode.</p> <p>A password is required to open the configuration utility.</p>
<p><b>DOC Service</b></p>	<p>Click <b>Enter</b> to perform DOC maintenance, including cleaning or replacing the Optical Adapter.</p>
<p><b>Field Service Tool</b></p>	<p>Click <b>Launch</b> to start the <b>Field Service Utility</b>.</p> <p><b>NOTE:</b> The service utility is for use by trained Field Service Engineers.</p> <p>A password is required to open the service utility.</p>
<p><b>External Monitor</b></p>	<p>This feature is not used on the OPTIS Integrated system.</p>

Table 10-16: Setup - Service Menu Settings (*continued*)

<b>Calibration Options</b>	The <b>Hide Calibration</b> checkbox controls the display of the calibration sequence when a catheter is first connected to the DOC. <b>Checked</b> - The calibration sequence is hidden. <b>Unchecked</b> - The calibration sequence is displayed on screen.
----------------------------	---

## Setup - Service/System Diagnostics Menu

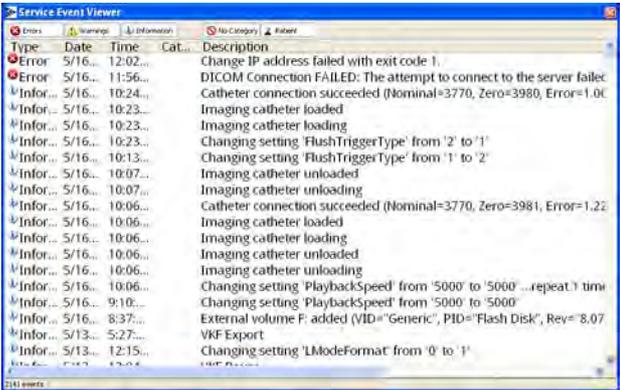The **System Diagnostics** menu in the **Setup** dialog box is provided to help St. Jude Medical Service monitor and analyze signal levels in the imaging engine either in real time or at the time a recording was made. When the tab is selected the system immediately displays the detected signal levels while in Acquisition mode or the signal values when the current image was captured when in Playback mode. If you encounter problems while using the system, you may be asked to access and report information from this tab

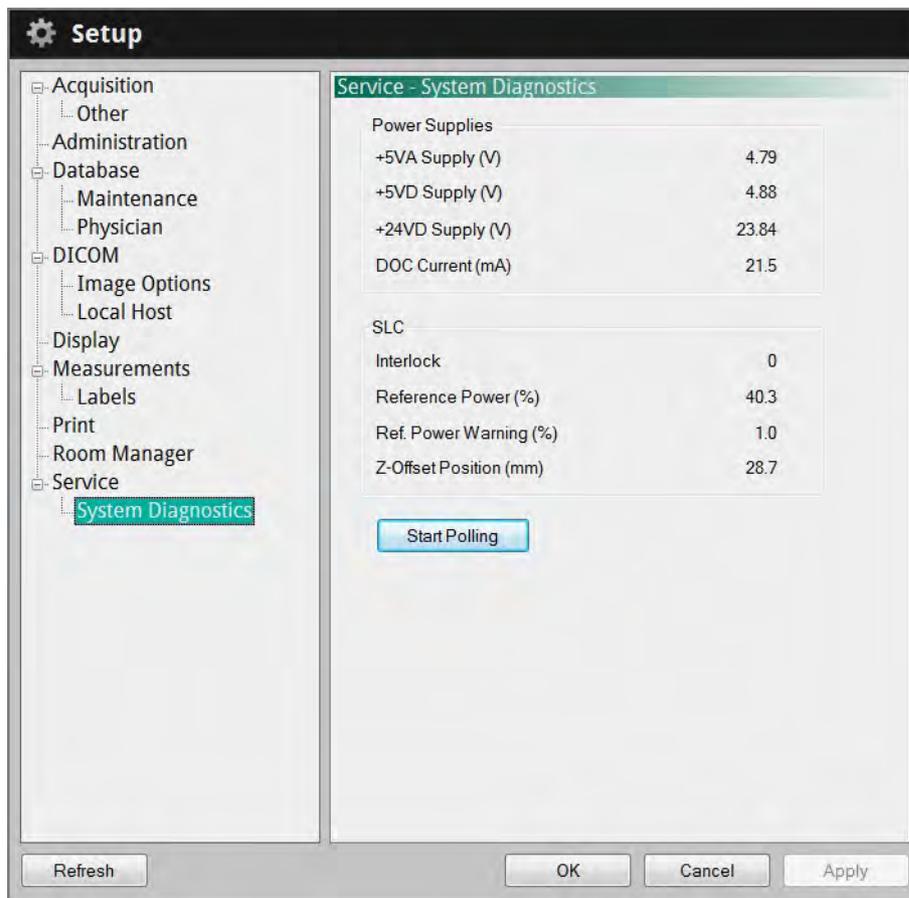


Figure 10-16: Setup - Service/System Diagnostics Menu

Table 10-17: Setup - Service/System Diagnostics Menu Settings

<b>Power Supplies</b>	<p><b>+5VA Supply</b> - The +5 volt AC power supply voltage.</p> <p><b>+5VD Supply</b> - The +5 volt DC power supply voltage.</p> <p><b>+24VD Supply</b> - The +24 volt DC power supply voltage.</p> <p><b>DOC Current</b> - The current being drawn by the DOC in mA.</p>
-----------------------	--

**User Interface Reference**  
**Setup - Service/System Diagnostics Menu**

---

Table 10-17: Setup - Service/System Diagnostics Menu Settings (*continued*)

<b>SLC</b>	<b>Interlock</b> - Displays the SLC interlock state. <b>Reference Power</b> - Displays the reference power value in %. <b>Ref. Power Warning</b> - Reference Power Warning value in %. <b>Z-Offset Position</b> - Z-Offset Position value in millimeters.
<b>Start Polling / Stop Polling</b>	Starts/stops automatic update of these values, in real time, every 50 milliseconds. <b>NOTE:</b> During playback, the <b>Polling</b> button is disabled, and the values represent the signal levels of the imaging engine at the image recording time.

Although the ILUMIEN OPTIS System conforms to laser emission standards and both international and European safety and electromagnetic compatibility standards, the system is intended for use only by medical personnel who have received ILUMIEN OPTIS System training. Only a trained operator can determine if ILUMIEN OPTIS System use is appropriate. An awareness of the system's limitations is essential to making that determination and assuring safe operation for both operator and patient.

This chapter includes:

- Precautions to assure patient and operator safety.
- How to avoid optical, electrical, explosion, and defibrillator hazards.
- How to make proper electrical connections.
- System imaging limitations.
- How to trace the source of electromagnetic interference.

**CAUTION:** Before using the ILUMIEN OPTIS System for the first time, be sure to read and understand all of the information in this chapter.

**NOTE:** The ILUMIEN OPTIS System complies with FDA performance standards for laser products except for deviations pursuant to laser Notice No. 50, dated July 26, 2001.

**NOTE:** LightLab Imaging, Inc. hereby declares that OPTIS Integrated system is in compliance with the essential requirements and other relevant provisions of Directive 1999/5/EC. A copy of the Full Declaration of Conformity can be obtained by contacting the EU Representative.

## **Patient Safety**

The ILUMIEN OPTIS System is intended for use only by medical personnel trained in its operation and skilled in the clinical procedures to be used.

To avoid any potential hazard to patients, follow the precautions outlined in this section.

**CAUTION:** Use only the Dragonfly Imaging catheters with the ILUMIEN OPTIS System. Use of other types of catheters may result in unsafe conditions for the patient and damage the ILUMIEN OPTIS System.

### **General**

**WARNING:** Failure to follow the guidelines described in these Instructions for Use and in the Instructions for Use provided with the accessories may result in injury to patients and damage to equipment.

- Only use Dragonfly™ Imaging Catheters. Always use under appropriate imaging guidance (endoscopy, x-ray fluoroscopy, or other appropriate guidance method).
- Only use PressureWire pressure transducers to report aortic distal pressure. Always use under appropriate imaging guidance (endoscopy, x-ray fluoroscopy, or other appropriate guidance method).
- Always read and follow the *Instructions for Use* supplied with the Dragonfly Imaging Catheter and with the PressureWire catheter.
- Always use controls, make adjustments and perform procedures as specified in these Instructions for Use.

### **Techniques to Minimize Patient Exposure**

The ILUMIEN OPTIS System meets the performance standards of laser-emitting products as established by IEC 60825-1. Although no harmful effects have been demonstrated for the near-infrared light wavelengths, intensities, and exposure times used during examinations with the ILUMIEN OPTIS System, St. Jude Medical recommends that you carefully read the warning labels on the system (see [Figure 11-1](#)) and follow these examination guidelines:

- Use OCT only when there is a good reason to do so.
- Use techniques that enable quick collection of clinical data and shorten procedure time.

## Operator Safety

### Avoiding Operator Light Emission Hazards

To avoid any potential light emission hazards to yourself or patients, adhere to the information provided in the safety labels that are located on the system (see [Figure 11-1](#)), and observe the precautions outlined in this section.

**WARNING: Failure to follow any of these precautions may cause possible serious damage to your eyes.**



Figure 11-1: Connector Panel Laser Safety Labels

- Avoid eye exposure. Do not look at or stare directly into the beam. Doing so may damage your eyes.
- Never view the laser output with optical instruments (for example, eye loupes, magnifiers, and microscopes). Doing so may damage your eyes.
- Use controls, make adjustments, and perform procedures only as specified in these Instructions for Use.

### Repetitive Strain Injury (RSI)

Repetitive use of a mouse and keyboard has been associated with Carpal Tunnel Syndrome (CTS) and related musculoskeletal problems. Follow these suggestions to help prevent these problems:

- Maintain your joints in optimum positions with a balanced posture, avoiding:
  - Static postures.
  - Exertion of force during repetitive motions.
  - Wrist flexion or deviation.
- Position the keyboard and monitor to minimize reaching and stretching.
- Take frequent breaks to give tissues time to recuperate from awkward positions and repetitive movements.

## **Moving the System**

When moving the system, observe these precautions:

**WARNING: Failure to follow any of these precautions may lead to a system tipping hazard, causing possible injury to people and damage to the system.**

- Be sure to turn off and disconnect the system cord from the wall outlet before beginning a move.
- Position the DOC cable appropriately before moving the system.
- The system weighs up to 95 kg (209 lbs) with all accessories installed, two people are required when moving it.
- Make sure that the system's wheels roll freely before beginning the move. Resolve any wheel problems before you move the system.

**NOTE:** Be sure the system brakes are in the up position (unlocked).

- To eliminate the potential danger of the system's tipping over, avoid ramps that are steeper than 5 degrees.
- Do not push the system by the monitors or monitor support mount.  A red circle with a diagonal slash over a black silhouette of a person pushing a cart.
- If you must move the system up or down ramps with an incline of more than 5 degrees, use two people.
- Do NOT lift a cart bearing the system to move it over uneven elevator entrances or other steps and barriers. Instead, find a route that avoids such problems.

**NOTE:** Wheelchair ramps usually have an incline of less than 5 degrees.

- When using a transport vehicle, be sure that it can handle the weight of the system components plus passengers.
- If a lift is used be sure the load capacity of the lift can accommodate the weight of the system components plus passengers.

## **Avoiding Electrical Hazards**

The isolation transformer in the ILUMIEN OPTIS System provides electrically isolated power to components supplied with the system or specified as part of the system.

**WARNING: All system components except the isolation transformer itself MUST be powered by, and only by, the isolation transformer in the ILUMIEN OPTIS System.**

This electrical isolation separates the system components, and thus the patient, from dangerous leakage currents. If any of the system's components are directly connected to a wall outlet or some power source other than the isolation transformer, the patient and the operator are no longer safely isolated and may be exposed to dangerous electric currents.

**WARNING: Failure to abide by the precautions detailed below causes the system and its use to be out of compliance with regulations and places the patient and the operator at risk of injury or death.**

- Do NOT use additional cables, extension cords or outlets with the ILUMIEN OPTIS System.

**NOTE:** The ILUMIEN OPTIS System's isolation transformer creates an electrically isolated Patient Interface for invasive imaging. This electrically isolated condition must be maintained throughout the procedure or during any contact with the patient.

- When transferring files inside the catheterization lab, use only port-powered USB drives.

**WARNING: Inside the catheterization lab only port-powered USB drives may be connected to the USB port. Connecting externally powered devices to the USB port in the patient vicinity may compromise electrical isolation and cause patient injury.**

**NOTE:** Outside the catheterization lab, IEC 60950-compliant, externally powered USB hard drives may be connected to the USB port.

- Do NOT remove system covers. Only qualified personnel should service the system. Accidentally contacting the electrical circuits inside the housing could cause serious injury.

## Making Proper Electrical Connections

Ensure the electrical connection for the system is properly rated (see [Figure 11-2](#)). Carefully follow the safety guidelines described in this section when connecting your system's power cord to the hospital or lab's AC outlet.

**WARNING:** Failure to follow the electrical connection precautions detailed in this section causes the system and its use to be out of compliance with regulations and places the patient and the operator at risk of injury or death and may damage the equipment.



Figure 11-2: Electrical Label

- Connect the system only to properly grounded (three-hole) hospital-grade AC outlets:
  - The circuit must accommodate an additional load of up to 400 VA.
- Replacement fuses are available through St. Jude Medical, part number: 11559-02 FUSE
- The power cord is to be used for disconnection from main power.
- Make sure that any devices that connect to the network interface of the ILUMIEN OPTIS System comply with the appropriate IEC/national standard and are certified to IEC 60950.
- Use no electrical peripherals within six feet of a patient unless the peripherals receive power from an isolation transformer that meets medical safety standards.

**NOTE:** If the ILUMIEN OPTIS System is used with peripherals that are powered from a separate wall outlet, then the combination is considered to be a Medical System. It is the user's responsibility to comply with IEC 60601-1-1 and test the Medical System according to the requirements.

## Explosion Hazard

**WARNING: Do NOT operate the ILUMIEN OPTIS System in the presence of flammable anesthetics. Doing so could lead to an explosion.**

## System Imaging Limitations

The ILUMIEN OPTIS System is intended for use by medical personnel who have received training in the use of the system. To determine if system use is appropriate, the trained user must be aware of system imaging limitations.

**CAUTION: Use only the Dragonfly™ Imaging catheters with the ILUMIEN OPTIS System. Other types of catheters may break if used.**

This section includes information about system capabilities and limitations for both vessel and tissue imaging.

### Considerations for Optimal Vessel Imaging

The ILUMIEN OPTIS System can be used to image through vessels or to image the inner surface of certain vessels. Since the Dragonfly catheter is smaller than the diameter of the vessel being imaged, the position of the catheter in the vessel has an effect on the portion of the vessel that can be imaged.

- Imaging range is greatest when the imaging catheter is centered in the lumen.
- Imaging range is least when the imaging catheter is placed eccentrically (off-center) in the lumen, against the wall of the vessel.

### Considerations for Optimal Tissue Imaging

The maximum imaging depth within a vessel wall ranges from approximately 0.9 mm to 1.3 mm and is limited by optical attenuation caused by scattering of the optical beam by microstructures in the vessel. For example, penetration of the OCT beam is deepest in calcified tissue and shallowest in dense fibrotic tissue.

## **Electromagnetic Compatibility**

The ILUMIEN OPTIS System is designed to meet the following electromagnetic compatibility standards:

- IEC 60601-1-2 (US and International).
- EN 60601-1-2 (Europe).

## **Electromagnetic Interference**

The system produce images by using digital signal processing techniques that operate in the radio frequency (RF) energy range. The system is therefore susceptible to interference generated by other RF energy sources such as medical devices, information technology products, or radio/television transmission towers. Tracing the source of radiated interference can be difficult.

In accordance with the standards identified in these Instructions for Use, no interference was observed. However, the trained user must determine if an artifact caused by radiated interference will negatively impact image quality and the subsequent study results.

To help identify the source of electromagnetic interference, ask the following questions:

- Is the interference intermittent or constant?
- Does the interference occur with one catheter only, or with other imaging catheters?
- Is the interference present if the system is moved to a different location in the facility? Examples: Placing the system close to an ECG cable can increase electromagnetic interference. Moving the cable or other medical equipment away from the system can reduce electromagnetic interference.

Please answer these questions before contacting your service representative. The answers will help a service representative determine if the problem is in the system or in the imaging environment.

## **Safety Functions Built Into the ILUMIEN OPTIS System**

The following safety functions have been built into the system:

- The system disables light output and disables all motors in these situations:
  - The optical fiber stops rotating due to mechanical failure.
  - Communication is lost between the imaging engine and the DOC.
  - Communication is lost between the computer and the imaging engine.
- Pressing the **Stop** button on the DOC disables power to the DOC and laser output (see [“The Drive-motor and Optical Controller \(DOC\)” on page 1-10](#)).

**Safety Information**  
**Safety Functions Built Into the Illumien Optis System**

---

## System - Safety & Regulatory

Table 12-1: System Safety & Regulatory Specifications

Category	Specifications
Regulatory Approvals	US 510(k) clearance. 
Safety standards system meets:	IEC 60601-1:2005 + CORR. 1 (2006) + CORR. 2 (2007)  EN60601-1: 2006, “Medical Electrical Equipment, Part 1: General Requirements for Safety”  EN60601-1-2: 2007, Electromagnetic Radiated Emissions Requirements for Medical Electrical Equipment - Group 1 Equipment, Class B for Non-Life Supporting Equipment  UL60601-1: 2003, “Medical Electrical Equipment, Part 1: General Requirements for Safety”  CAN/CSA C22.2 No. 60601-1:2008, “Medical Electrical Equipment, Part 1: General Requirements for Safety”  IEC 60825-1, 2nd Ed., 2007: Safety of Laser Products
Electromagnetic compatibility (EMC)	Refer to <a href="#">Table 12-5</a> , <a href="#">Table 12-6</a> , <a href="#">Table 12-7</a> , and <a href="#">Table 12-8</a> for detailed specifications.

## System Specifications

### System - Safety & Regulatory

Table 12-1: System Safety & Regulatory Specifications (*continued*)

Category	Specifications
<b>Classifications</b>	
Type of protection, shock	Class 1
Degree of protection, shock	Type CF <ul style="list-style-type: none"> <li>• DOC with catheter (CF label at DOC cable exit on Connector Panel)</li> </ul>
Degree of protection, ingress	<b>Console</b> - IPX0 <b>DOC</b> - IPX0, use with Sterile DOC cover for ingress protection
Method of Disinfection	Console and DOC will withstand without damage or deterioration disinfection by wiping with common hospital disinfectants including Cidex (Glutaraldehyde 3.4%).
Flammable mixtures	Not for use in presence of flammable anesthetic mixture with air or with oxygen or nitrous oxide.
Mode of operation	Continuous
<b>Leakage &amp; Auxiliary Current</b>	
Chassis leakage current	< 100 $\mu$ a rms normal condition < 500 $\mu$ a rms single-fault condition
Patient leakage current	Measured at patient end of DOC: < 10 $\mu$ a rms normal condition < 50 $\mu$ a rms single-fault condition
<b>Software Safety Features</b>	
The computer and software are designed with the following security features. These features do not require any user configuration or action.	
<ul style="list-style-type: none"> <li>• Prevents unauthorized access to the operating system.</li> <li>• Prevents installation or execution of unauthorized software</li> <li>• Prevents infiltration via Ethernet connection.</li> <li>• Performs data integrity check when archiving data to external media.</li> </ul>	

## System - Electrical and Physical

Table 12-2: System Electrical and Physical Specifications

<b>Parameter</b>	<b>Specification</b>
<b>Power Input</b>	
Line voltage	100/120/220/240 VAC $\pm$ 10%, user selectable 50/60 Hz $\pm$ 1 Hz
Power consumption	Active: < 400 VA Standby: < 30 VA
<b>Transport and Storage Conditions (Permissible ranges)</b>	
Ambient temp	-25 to +50 degrees C
Relative humidity	10% - 95%, including condensing
Atmospheric pressure	500 to 1060 mBar
<b>Operating Conditions</b>	
Ambient temperature	+10 to +32 degrees C
Relative humidity	10% to 85%, non-condensing
Atmospheric pressure	700 to 1060 mBar
<b>Mechanical Specifications</b>	
Weight	95 kg (209 lbs) max with all accessories
Overall Dimensions	145 cm H x 61 cm W x 71 cm D $\pm$ 5 mm

## Imaging Specifications

Table 12-3: Imaging Specifications

<b>Parameter</b>	<b>Specification</b>
<b>Optical Parameters - Measured at System Aperture (DOC Optical Port)</b>	
Scanning Laser Source Optical Power	22.6 mW maximum @ 1305 nm $\pm$ 55 nm (Class 1M Laser Output per IEC 60825-1)
Visible Laser Optical Power	1.45 mW maximum @ 670 nm (nominal) (Class 1M Laser Output per IEC 60825-1)
<b>Pullback Parameters</b>	
Pullback Range	75 mm  (If connected to a C7 Dragonfly catheter, the range is 54 mm.)
Pullback Speed Settings	18.0 mm/sec, 36.0 mm/sec  (If connected to a C7 Dragonfly catheter, the speed settings are 10.0 mm/sec, 20.0 mm/sec, 25.0 mm/sec.)
<b>General Scan Parameters</b>	
A-Scan Range in Air	7.0 mm
A-Scan Range in Contrast	4.83 mm
Diameter Measurement Accuracy	7% $\pm$ 0.1 mm
Area Measurement Accuracy	10% $\pm$ 0.1 mm <sup>2</sup>
Axial Resolution	$\leq$ 20 $\mu$ m in tissue
Optical Sensitivity	90 db minimum
A-Scans per second	90 kHz (nominal)
Frame Rate	180 frames/second (Hz)  (If connected to a C7 Dragonfly catheter, the Frame Rate is 100 frames/second (Hz).)

## FFR Specifications

Table 12-4: FFR Specifications

<b>Parameter</b>	<b>Specification</b>
<b>AO Pressure (Wi-Box to ILUMIEN OPTIS System)</b>	
Operating pressure	-200 to +450 mm Hg
Accuracy	+/- 1 mm Hg or +/- 1% of reading, whichever is greater
<b>PW Pressure</b>	
Operating pressure	-30 to +300 mm Hg
Accuracy	$\pm 1$ mm Hg plus $\pm 1\%$ of reading (-30 to 50 mm Hg) $\pm 3\%$ of reading (50 to 300 mm Hg)
<b>AO Pressure (Wi-Box to hemodynamic recording system)</b>	
Direct galvanic connection	
Max pressure shift	<2 mm Hg
<b>Radio Specification</b>	
Frequency range	2.4000 - 2.4835 GHz
Type	Frequency hopping spread spectrum (FHSS)
Range	0 - 4 m
Delay time	<20 ms

## **Electromagnetic Emissions**

Table 12-5: Guidance and Manufacturer's Declaration - Electromagnetic Emissions

The ILUMIEN OPTIS System is intended for use in the electromagnetic environment specified below. The customer or user of the ILUMIEN OPTIS System should assure that it is used in such an environment.

<b>Emissions Test</b>	<b>Compliance</b>	<b>Electromagnetic Environment - Guidance</b>
RF emissions CISPR 11	Group 1	The ILUMIEN OPTIS System uses RF energy only for its internal function. Therefore, its RF emissions are very low and not likely to cause any interference with nearby electronic equipment.
RF emissions CISPR 11	Class B	The ILUMIEN OPTIS System is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic proposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	

## Electromagnetic Immunity

Table 12-6: Guidance and Manufacturer's Declaration - Electromagnetic Immunity

The ILUMIEN OPTIS System is intended for use in the electromagnetic environment specified below. The customer or user of the ILUMIEN OPTIS System should assure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	± 2 kV for power supply lines ± 1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV lines(s) to line(s) ± 2 kV lines(s) to earth	± 1 kV lines(s) to line(s) ± 2 kV lines(s) to earth	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	< 5 % $U_T$ (>95% dip in $U_T$ ) for 0.5 cycle 40% $U_T$ (>60% dip in $U_T$ ) for 5 cycles 70% $U_T$ (30% dip in $U_T$ ) for 25 cycles <5% $U_T$ (>95% dip in $U_T$ ) for 5 seconds	< 5% $U_T$ (>95% dip in $U_T$ ) for 0.5 cycle 40% $U_T$ (>60% dip in $U_T$ ) for 5 cycles <70% $U_T$ (>30% dip in $U_T$ ) for 25 cycles <5% $U_T$ (>95% dip in $U_T$ ) for 5 seconds	Mains power quality should be that of a typical commercial or hospital environment. If the user of the ILUMIEN OPTIS System requires continued operation during power mains interruptions, it is recommended that the ILUMIEN OPTIS be powered from an uninterruptible power supply or battery.

**NOTE:**  $U_T$  is the AC mains voltage prior to application of the test level.

## System Specifications

### Electromagnetic Immunity

---

Table 12-6: Guidance and Manufacturer's Declaration - Electromagnetic Immunity

The ILUMIEN OPTIS System is intended for use in the electromagnetic environment specified below. The customer or user of the ILUMIEN OPTIS System should assure that it is used in such an environment. (*continued*)

---

<b>Immunity Test</b>	<b>IEC 60601 Test Level</b>	<b>Compliance Level</b>	<b>Electromagnetic Environment - Guidance</b>
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical commercial or hospital environment.

---

Table 12-7: Guidance and Manufacturer’s Declaration - Electromagnetic Immunity

The ILUMIEN OPTIS System is intended for use in the electromagnetic environment specified in this section. The customer or user of the ILUMIEN OPTIS System should assure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
			Portable and mobile RF communications equipment should be used no closer to any part of the ILUMIEN OPTIS System, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
Conducted RF IEC 61000-4-6	3 Vrms  150 kHz to 80 MHz	3 Vrms	Recommended separation distance $d = 1.2\sqrt{P}$ $d = 1.2\sqrt{P}$ 80 MHz to 800 MHz $d = 2.3\sqrt{P}$ 800 MHz to 2.5 GHz
Radiated RF IEC 61000-4-3	3 V/m  80 MHz to 2.5 GHz	3 V/m	Where P is the maximum output power rating of the transmitter in Watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).  Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey <sup>a</sup> , should be less than the compliance level in each frequency range <sup>b</sup> .  Interference may occur in the vicinity of equipment marked with this symbol:  

**NOTE 1:** At 80 MHz and 800 MHz, the higher frequency range applies.

**NOTE 2:** These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

## **System Specifications**

### **Electromagnetic Immunity**

---

a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the ILUMIEN OPTIS System is used exceeds the applicable RF compliance level above, the ILUMIEN OPTIS System should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the ILUMIEN OPTIS System.

b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

## Recommended Separation Distances

Table 12-8: Recommended separation distances between portable and mobile RF communications equipment and the ILUMIEN OPTIS System

The ILUMIEN OPTIS System is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the ILUMIEN OPTIS System can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the ILUMIEN OPTIS System as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1.2\sqrt{P}$	80 MHz to 800 MHz $d = 1.2\sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3\sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3

For transmitters rated at a maximum output power not listed above, the recommended separation distance  $d$  in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where  $P$  is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

**NOTE 1:** At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

**NOTE 2:** These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

## **FCC Statement**

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.

Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by St. Jude Medical could void the user's authority to operate the equipment.

## **Essential Performance Defined by Operating Mode**

| During testing, the OCT Imaging System (MN: C8/OPTIS Integrated) was operating as follows: SSOCT Data Viewer is displaying white noise. The CardioUI Development software is running. The probe is active and displaying an image. Image on screen is stable and without error. Also host PC is pinging the EUT.

# Index

---

## Symbols

---

%AS, [7-2](#), [7-10](#)

%DS, [7-2](#), [7-12](#)

## Numerics

---

3D Display, [7-23](#)

## A

---

Acquisition mode, [10-29](#)

Adding Points, [7-18](#)

Annotation

adding text, [7-8](#)

text, [7-9](#)

Area Measurements

automatic, [7-6](#)

Auto Cycle Drawing Color, [10-23](#)

Automatic MLA and %DS option, [7-20](#)

## B

---

Black Level, [9-15](#)

Bookmark

clear all bookmarks, [6-10](#)

creating, [6-10](#)

next bookmark, [6-10](#)

previous bookmark, [6-10](#)

## C

---

Calculation

%AS, [7-10](#)

%DS, [7-12](#)

Callout

adding, [7-9](#)

position, [7-9](#)

Callouts

adding, [7-8](#)

Catheter

connection to DOC, [4-5](#), [5-8](#)

disconnecting from DOC, [5-18](#)

insertion, [5-13](#)

positioning, [5-13](#)

preparation, [5-7](#)

purge, [5-9](#)

removal, [5-18](#)

stop movement, [5-19](#)

Caution, meaning and format, [Front-iv](#)

Cleaning, [9-2](#)

Complications, [1-16](#)

Complications from Use, [1-16](#)

Connection

network, [1-8](#)

power, [2-2](#)

USB, [1-8](#)

video, [1-8](#)

Connector Panel

location, [1-4](#)

Contact Information

company, [Front-ii](#)

service, [9-2](#)

Contraindications for Use, [1-13](#)

Control Point Size, [10-23](#)

Conventions Used in Manual, [Front-iv](#)

Creating a New Patient, [3-6](#)

Cross-section View

zoom, [7-15](#)

## Index

### D

---

### D

---

#### Database

- create a new patient, [3-6](#)
- import, [8-20](#)
- Setup dialog box, [10-7](#)
- statistics, [8-26](#), [10-9](#)

#### Database Tab, [10-7](#)

#### Delete

- files, [8-23](#)
- measurements, [7-19](#)
- points, [7-18](#)

#### Depth Calibration Marks, [5-16](#)

#### Diagnostics Tab, [10-29](#)

#### DICOM

- viewer, [8-16](#)

#### DICOM Tab, [10-12](#)

#### DOC

- catheter connection, [4-5](#), [5-8](#)
- cleaning optical connection, [9-5](#)
- description, [1-10](#)
- preparation, [5-6](#)
- replacing optical adapter, [9-8](#)

#### Drive-motor and Optical Controller, [1-10](#)

### E

---

#### Electrical Connections, [11-6](#)

#### Electrical Hazards, [11-5](#)

#### Enable button on DOC, [1-10](#)

#### External Drive, [8-9](#), [8-10](#), [8-12](#)

### F

---

#### FFR Procedure

- catheter connection, [4-5](#)
- materials and equipment, [4-1](#)

#### File Size, [8-6](#)

#### Files

- deleting, [8-23](#)

#### Flush Medium

- warning, [1-14](#)

### G

---

#### Green's Theorem, [7-6](#)

### H

---

#### Hazards

- electrical, [11-5](#)
- explosion, [11-7](#)
- light emission, [11-3](#)
- repetitive strain injury, [11-3](#)

#### Home menu, [3-2](#)

### I

---

#### Image Compression

- caution, [7-1](#), [8-7](#)

#### Image Format

- raw, [8-4](#)
- standard, [8-5](#)

#### Imaging

- limitations, [11-7](#)
- tissue, [11-7](#)
- vessel, [11-7](#)

#### Imaging Specifications

- optical parameters, [12-4](#)
- pullback parameters, [12-4](#)
- scan parameters, [12-4](#)

#### Import OCT Files, [8-20](#)

#### Indications for Use, [1-12](#)

#### Infection Control, [9-13](#)

#### Instructions for Use

- other manuals, [Front-iv](#)

#### Intended Use, [1-12](#)

### L

---

#### Length Measurement, [7-5](#)

#### Line Width, [10-23](#)

#### Live Mode, [10-29](#)

#### Live View, [5-3](#)

#### Live View Button on DOC, [1-10](#)

---

**L-Mode**

- caution, 7-3
- cut-plane, 5-16
- limitations, 6-5
- measurements and annotations, 7-3

Lock LED on DOC, 1-10

Log Files, 9-11, 10-26

Lumen Profile, 7-20

---

**M**

Maintenance, 9-4

Manual

- conventions, *Front-iv*

Measurement Accuracy, 7-3

Measurement and Annotation tools, 7-2

Measurements

- adding points, 7-18
- caution, 7-1, 8-7
- deleting all measurements, 7-19
- deleting individual measurements, 7-19
- deleting measurements, 7-19
- deleting points, 7-18
- editing, 7-17
- length, 7-5
- moving individual points, 7-18

Minimum Lumen Area, 7-21

MLA, 7-21

Monitor Setup, 2-6

Moving Points, 7-18

Moving System, 11-4

---

**N**

Near-infrared Light, 1-1, 11-2

Network Connection, 1-8

Note, meaning and format, *Front-iv*

---

**O**

OCT Database, see *Database*

**OCT Procedure**

- catheter connection, 5-8
- catheter insertion, 5-13
- catheter positioning, 5-13
- catheter preparation, 5-7
- completing procedure, 5-18
- DOC preparation, 5-6
- materials and equipment, 5-1
- recording, 5-16

**Operator Safety**

- light emission hazards, 11-3
- moving the system, 11-4
- repetitive strain injury, 11-3

Optical Coherence Tomography, 1-1

Optical Parameters, 12-4

Optimal Tissue Imaging, 11-7

Optimal Vessel Imaging, 11-7

---

**P****Patient**

- acquiring image, 5-16
- creating a record, 3-6
- minimizing exposure, 11-2
- safety, 11-2

**Patient Entry**

- creating, 3-6

Patient Record

- Open, Create, 3-1

Pen Color, 10-23

Percent Area Stenosis, 7-2, 7-10

Percent Diameter Stenosis, 7-2, 7-12

PIU. See *DOC*

Playback Calibration

- setting, 6-7

Playback Mode, 10-29

playback range, 6-11

Positioning the System, 2-1

Power

- off, 2-4
- on, 2-3

Power In, 12-3

Precautions for Use, 1-15

## Index

### R

---

#### Pullback

- parameters, [12-4](#)
- stop, [5-19](#)
- trigger, [10-3](#)

#### Pullback motion LEDs on DOC, [1-10](#)

#### Purge Catheter, [5-9](#)

### R

---

#### Range

- playback, [6-11](#)

#### Raw Format

- description, [8-4](#)

#### Recording

- status, [5-16](#)

#### Recording calibration marks, [5-16](#)

#### Remote DICOM Store, [8-11](#)

#### Reviewing

- saved images, [6-1](#)

### S

---

#### Safety

- functions, [11-9](#)
- operator, [11-3](#)
- patient, [11-2](#)

#### Segmented Lumen, [7-25](#)

#### Select Measurement

- percent area stenosis, [7-10](#)
- percent diameter stenosis, [7-12](#)

#### Select Patient menu, [3-2](#)

#### Service Tab, [10-26](#)

#### Setup Dialog Box

- Database tab, [10-7](#)
- Diagnostics tab, [10-29](#)
- DICOM tab, [10-12](#)
- Service tab, [10-26](#)

#### Shutdown

- procedure, [2-4](#)

#### Specifications

- electrical and physical, [12-3](#)
- electromagnetic, [12-6](#), [12-7](#), [12-9](#)
- imaging, [12-4](#)
- pressure and radio, [12-5](#)
- safety and regulatory, [12-1](#)

#### Standard File Format, [8-6](#)

#### Standard Format

- description, [8-5](#)

#### Standby View, [5-3](#)

#### Stop Button on DOC, [1-10](#)

#### System

- components, [1-2](#)
- features, [1-1](#)
- imaging limitations, [11-7](#)
- moving, [11-4](#)
- power on, [2-3](#)
- setup, [2-1](#)
- shut down, [2-4](#)

#### System Display

- during image acquisition, [5-16](#)

#### System Specifications

- classifications, [12-2](#)
- leakage & auxiliary current, [12-2](#)
- mechanical specifications, [12-3](#)
- operating conditions, [12-3](#)
- power input, [12-3](#)
- software safety features, [12-2](#)
- transport and storage conditions, [12-3](#)

### T

---

#### Text

- adding, [7-8](#)
- display options, [7-9](#)
- entry, [7-9](#)
- position, [7-9](#)

#### Tissue Imaging, [11-7](#)

#### Transfer Messages, [8-24](#), [8-26](#)

#### Trigger Type, [10-3](#)

#### Troubleshooting

- connections, [9-15](#)
- DOC, [9-15](#)
- general, [9-14](#)
- imaging, [9-15](#)
- PressureWire
  - power failure, [4-15](#)

### U

---

#### Unload button on DOC, [1-10](#)

## USB

- connection, [1-8](#)
- media, [8-2](#)

## V

---

- Vessel Imaging, [11-7](#)
- Video Connection, [1-8](#)

## W

---

- Warning, [1-14](#)
- Warning, meaning and format, [Front-iv](#)
- Weight, [12-3](#)
- White Level, [9-15](#)
- Wi-Box
  - cathlab installation, [1-11](#)

## Z

---

### Zoom

- DICOM Viewer, [8-17](#)
- measurements, [7-3](#)
- region, [7-15](#)





