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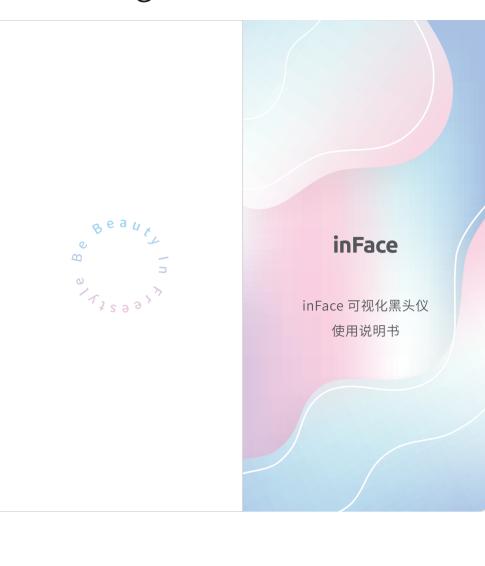
园U6栋807 官方网站: www.migoosmart.com

Welcome to use inFace

inFace Visual Blackhead Remover is independently developed and manufactured by Shenzhen inFace Technology Co., Ltd., which is affiliated to Shenzhen Migoosmart Co., Ltd.. create products with technology and intellectual, and striving

Precautions ----- 12 FCC Warnning ----- 13 Product Details ------ 14 Product Parameters ----- 15 Instructions ----- 15 Operations ------ 16 Charging Instructions ----- 18 Cleaning and Maintenance —----- 18 Q&A ----- 18 Environmental Protection ----- 19 Waste Treatment —----- 20

Product List

Funcitions of the Product —----- 12

the product and keep it properly.

After-sales Protection and Service —----- 20 Certification ------ 21 Warranty Card ----- 21 After-sales Contact Information ----- 21

is product out of the reach of children.

inFace Visual Blackhead Remover Functions of the Product

Removing Blackheads | Exfoliating Cleaning Pores | Blue Light Anti-inflammatory

Precautions Keep the product away from bathtub, shower basin, washba-sin or other water-containing utensils, and keep the product and the power adapter dry. Dry your hands to prevent water or other liquids from it may cause the skin to become red, purple and swollen. It is recommended that the detention time in the same area during use should not exceed 3 seconds, and each treatment time 6. Do not use this product on skin with open wounds, acne or 7.Do not use it on cuts, areas that have undergone cosmetic

eating phenomenon, it is not a product failure, please rest assured to use it.

10.After using the product, the skin may appear reddish condition, which is a normal phenomenon, generally within 1 hour the redness will self-fade.
11.Do not look directly at the running LED light at close range. 12. Do not allow children to use this product, please place

been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
-children shall not play with the appliance FCC Warnning

party responsible for compliance could void the user's authority to operate the equipment.

the FCC Rules. These limits are designed to provide reasonal guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference

Reorient or relocate the receiving antenna.

Any Changes or modifications not expressly approved by the

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including protection against harmful interference in a residential installation. This equipment generates uses and can radiate

to radio or television reception, which can be determined by to correct the interference by one or more of the following -Increase the separation between the equipment and receiver.
-Connect the equipment into an outlet on a circuit different -Consult the dealer or an experienced radio/TV technician for

This device complies with part 15 of the FCC Rules. Operation is


2.Suction Nozzle x4 3.Silicone O Ring x2 4.Dust Net x2 5.USB Charging Cable x1
6.User Manual and Warranty Card x1
(If one or more of the above accessories is missing, please

contact our after-sales customer service center.) ---

-this appliance can be used by children aged from 8 years and above and persons with reduced physical,sensory or mental capabilities or lack of experience and knowledge if they have

13. If the product is damaged, please stop using it.

14. If you have any skin problems, please consult a doctor before using this product.

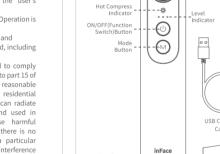inFace Visual User Manual and Blackhead Remover Warranty Card

requirement. The device can be used in portable exposure

Product Details

Model: CF-DSE Battery Capacity: 800mAh
Input Voltage: 5V = Input Power: 5W
Rated Voltage: 3.7V = Rated Power: 6W
Standby Current: About 30uA
Charging Time: About 2 hours
Charging Current: 600mA
Product Size: 160°31.5°31.5mm Net Weight: About 120g
Continuous Working Time: About 16 mins Working Temperature Range: 0°C~40°C 1.Long press the "ON / OFF(Function Switch) Button " for 1.5 seconds to turn on the device. Level 1-white LED flashes to to switch between "hot compress-blue light-cleaning" and

Product Parameters

Product Name: inFace Visual Blackheac Model: CF-05E Batter

the device will automatically shut down; Switch) Button " for 1.5 seconds, the device will shut down 5.In the "cleaning mode", short press the " ON / OFF(Function Switch) Button " to switch between "level 1 - level 2 - level 3 - standby", and the corresponding indicator lights are device is heated to 40 \pm 2 degrees and kept at a constant temperature, and the indicator light of "hot compress mode" 7.In the "blue light mode", the blue LED at the top of the device lights up; 8.When the battery voltage is lower than 3.4V, the low battery

mode. If there is no operation within 90 seconds of standby,

red indicator flashes for 90 seconds and automatically shuts

9.When charging, the three white indicators (LED1 / LED2 /

/ LED2 / LED3) are always on after fully charged; 10.Do not use when charging, the key operation is invalid.

3 4

[Steps for Usage] Download and install the APP (see APP download and install steps for details);

2.Connect the device with the APP (see APP function description for details); 3.Clean the face before using, apply a hot towel or use a the pores, and then apply a blackhead export nasal mask, turn on the device - hot compress function to accelerate the softening of the blackheads and help them easily suck out 4.Select the suction nozzle suitable for your needs, install it on the top of the device and test the suction force on the palm or forearm before using it; 5.When using, you can select suitable level according to your skin condition by pressing " ON / OFF(Function Switch) When using, make sure that the suction head is close to the

Use your mobile phone browser to the QR code in the manual, then

15 16

10.Facial cleaning care recommended once a week, the specific frequency of use must be adjusted according to their

Connecting to the camera of the device through APP makes it more intuitive and convenient to clean pores during use. 1. Turn on the device and enter the standby state, then click on

5 6

the phone settings page to select the WIFI name corresponding to the device (such as inFace-xxxx), and connect. 2.Open the " UDU " APP and go to the home page. If the device is not connected, the APP will prompt for operation instruc-

skin, and there should be no gaps, otherwise vacuum negative pressure will not be formed, and the blackheads will not be effectively sucked out;
7. Keep moving when using the device, and the detention time 8.After cleaning, you can switch to " blue light mode " for skin 9.After using, long press the "ON / OFF(Function Switch)
Button " for 1.5 seconds to turn off the device;

iOS system and Android system are common

Turn on 4.During using, the working position of the blackhead met can be adjusted appropriately according to the condition the displayed facial skin. 5.After clicking " Video ", the images transmitted by the device camera can be saved to the mobile phone as a video.
6.After clicking " Photo ", the image transmitted by the camera

Photo Button / Video Button

level based on the skin conditions.

7.After using it, just close the APP in the background of the

Charging Instructions L.Please use a standard qualified DC5V/1A adapter for

always on.

4.When the battery is low, the red indicator will flash slowly for 90 seconds. Please use it after fully charged.

** Connecting to the power supply, if charging is not started after 10 seconds, please confirm whether the Type-C cable is connected badly or whether the power supply is normal

corrosive solution wipe or wash with water. 4. Device camera lens with a cotton swab for cleaning, do not touch the lens with a hard object to avoid scratching the Lens.

5. Please store the device in the normal drying room.

3.After being fully charged, the three level indicators are

Cleaning and Maintenance When cleaning, make sure that the device is turned off and cannot be washed with water to avoid electric shock. When cleaning the suction nozzle, please remove the suction 3.The device is not waterproof, after each use can be used wet

an appropriate suction nozzle and select appropriate suction

Q: How long can it be using once it is fully charged? A: It can be using continuously for 60 mins when it is connected to the mobile phone and for 90 mins when it is not connect-

A: ① Due to the suction force of vacuum negative pressure, the skin will be slightly flushed following suction, which is normal and will fade away within 1 hour. ② It is recommended not to keep the device on one skin area for more than 3 sec when using the device. Ice compress can be conducted following

A: ① Please apply a lower level of suction; ② Please select a smaller suction nozzle; ③ Please keep moving during operation, Each treatment time should not exceed 5 minutes. Q: How often do you use it?

O: Why can't the APP be connected?

ate the recycling of useful materials.

A: Once in a week for sensitive skin and dry skin, twice a week for oily and mixed skin. Q: How to select an appropriate level? A: Before using the device, please check the suction force with your palm or forearm, then for the face, it is recommended to start the test from "Level 1", and then gradually adjust it to the

A: ① Please ensure that the visual blackhead vacuum device is turned on. ② WIFI in setting in your phone, then connect to the right WIFI of the device. Q: Why does the APP always about crash back? A: Probably due to low battery, please try again once the device is fully charged.

A: ① In normal cases, pores are closed. Before using the device, it is recommended to apply hot compress or use the facial steamer to open the pores, and then use the blackhead Environmental Protection 1.This product and packaging materials and accessories are recyclable materials.
 2.Divide waste into different recycling levels to facilit-

Waste Treatment

battery; the battery must be safely disposed.

end of its useful life. Please consult the local govern

After-sales Protection and Service

7 8

1.inFace provides a one-year warranty for products service standard. The warranty period starts from the date of 2.During the warranty period, for the failure of the product due to defects in the material or production process, we wil replace the parts or the entire product as appropriate. 3.To apply for warranty, the product and the proof (purchase invoice) that can prove the purchase date must be handed

19 20

outlet, man-made damage, etc.
②Product aging and wear caused by normal use. ries, resulting in product problems. 4 Consumables, accessories, etc..

Shenzhen Migoosmart Co., Ltd. After-sales Contact Information ①Damage caused by improper use, for example: using tl Tel: 0086-400-007-8388 wrong power supply, inserting the product into an unsuitable

⑤Damage caused by force majeure factors (such as lightning earthquake, flood, etc.). Note: Do not immerse the instrument in water at any time Damage to the machine caused by human factors i covered by the warranty. If there is any change in the above list, the relevant regulations newly formulated by the company or the relevant laws and regulations updated by the country shall prevail. For specific warranty terms, please consult the