


clarisonic

TOUCH

Sonic Infusion

USER GUIDE **DRAFT**



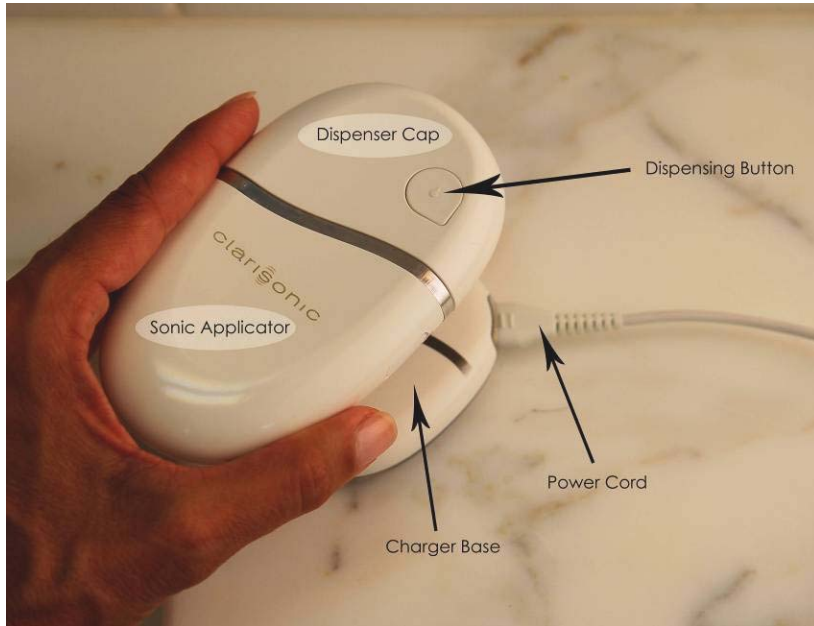

Anti-Aging Serum

Scientifically developed and optimized for the sonic applicator, this antioxidant-rich Sea Serum provides a multi-layered approach to balance and protect the skin from age-related and environmental damage. Nutrient and polysaccharide-rich marine and botanical extracts create a skin-firming complex designed to:

- Create a moisture protective complex for immediate and long-lasting hydration
- Promote cellular metabolism and collagen synthesis
- Reduce the appearance of fine lines and wrinkles
- Soothe and calm irritation-prone skin
- Brighten and even skin tone

Clarisonic Touch

Clarisonic Touch is a sonic applicator with a companion Sea Serum scientifically developed and optimized to help correct the skin's visible signs of aging, particularly targeting fine lines and wrinkles in the area around your eyes (commonly referred to as "crow's feet"). The applicator uses a gentle micro-pulsing action to effectively apply a measured amount of serum to the area(s) you need it most. The sonic action allows more effective absorption of the serum, immediately reducing the appearance of fine lines and wrinkles and promoting long-lasting, cumulative results.

The Sonic Touch

Tapping at more than 125 strokes/second, the sonic applicator works with the natural elasticity of your skin to replicate traditional "ring finger tapping" of anti-aging creams to fine lines and wrinkles. Touch's gentle kneading action allows the serum to reach the layers of the epidermis where it can work to reveal smooth, healthy, hydrated skin.

READ ALL INSTRUCTIONS BEFORE USE

Using the Touch

Prior to the first use, charge the Clarisonic Touch for 24 hours. See *Charging*.

! ALLERGY WARNING

Perform skin test at least 48 hours before using.

1. Apply single dose of serum to inside of forearm.
2. Rub excess serum onto arm.


If bumps, burning, skin discoloration, or other allergic reaction occurs, discontinue use and consult your physician.

3. Do not use on open sores or pustular acne.

Touch is gentle for use twice daily, morning and evening, on cleansed skin. Use as a skin-and-eye refresher when needed, even over facial makeup.

Cleanse your skin with the Clarisonic Skin Cleansing brush (or follow your manual cleansing ritual). After cleansing:

1. Remove the dispenser cap (containing the serum) from the Sonic Applicator.
2. Place the nozzle (located on the inside surface of the cap) slightly above the tip of the Touch.
3. Press and release the button on the side of the dispenser cap to dispense the serum.

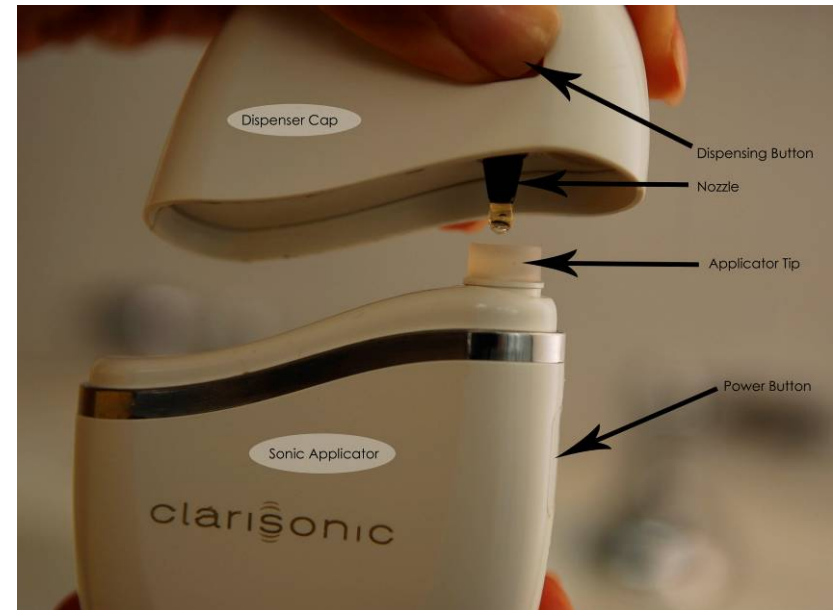

4. Place the applicator tip flat to the skin on the crow's feet area of the eye.
5. To activate the Touch, press and release the power button located on the side.
6. Follow the rim of the bone outlining the eye area to apply the serum using light pressure and small, overlapping, circular strokes.
 - ⇒ **Do not** use the Touch on the eye lids or soft tissues directly around your eye.
 - Touch will automatically turn off after 20 seconds.



8. Continue to use the Touch as outlined above to apply the serum to the other eye.
9. After use gently wipe the applicator tip with a soft, damp cloth, and replace the cap for storage.
 - ⇒ **Do not** immerse the Touch or place under running water.

NOTE: Any serum remaining on the treated areas may be gently tapped into the skin using your ring finger, if needed.

7. Repeat steps 2 through 6 to treat the next area.

Charging the Touch

! DANGER – Electrocution Risk with Charger

- ⇒ Keep charger away from water immersion.
- ⇒ If charger falls in water, do not touch it! Unplug from wall immediately.
- ⇒ Never use while bathing.
- ⇒ Do not use or store where charger can fall or be pulled into tub or sink.

Your Clarisonic Touch is rechargeable with the included charger, which magnetically attaches to the applicator. A fully charged Touch provides 30 days of use when used twice per day x two 20-second uses.

Charging Instructions

1. Attach the power cord to the charger.
2. Plug the cord into a working outlet.
3. Place the Touch onto the charger as shown below. (Logo facing up, applicator chrome band directly over charger gray band).
 - The button of the dispensing cap should align with the power cord

The Touch will magnetically attach to the charger.



Charge Indicators

Action	Indicator
Attach Touch to charger	Magnetic attachment. If the applicator slides off the charger, it is not properly attached.
Charge the applicator	Blue LED flashes while charging Blue LED glows solid when fully charged
After xx minutes of use, the applicator needs to be charged	[Define indicator] Allow 24 hours for a full charge.

The charger is compatible with 100-240 VAC; 50/60 Hz, allowing it to be used worldwide. If the shape of the plug does not fit the power outlet, use an attachment plug adapter of the proper configuration for the power outlet.

Caring for the Touch

To keep your Touch in optimum condition, care for it as follows:

- Do not use abrasive cleansers or chemicals to clean any part of the Touch.
- Do not put any part of the Touch in a dishwasher or immerse in water.

APPLICATOR TIP

After use gently wipe the applicator tip with a soft, damp cloth.

The applicator tip should be replaced after 90 days of use. Each new Serum Replenishment package includes a new applicator tip. To replace the tip, gently pop the tip from the applicator. Do not pull on the soft tip material.
[need illustration/photo of this]

APPLICATOR

The applicator is sealed. Never attempt to open the applicator for any reason. The applicator is water-resistant, but do not immerse in water.

REPLACING THE DISPENSING CAP

Each dispensing cap provides approximately 45 days of serum for sonic infusion when used twice/day in the eye area. When the dispensing cap is nearing empty, you may have to press the button more than once to dispense serum.

The dispensing cap cannot be refilled. Please recycle empty dispensing caps through our recycle program. The Serum Replenishment package includes a pre-paid return envelope for your convenience in recycling dispensing caps.

Storing the Touch

Clarisonic Touch can be safely stored between 32° and 120°F. Subjecting Touch to extreme heat or cold could result in damaging the unit. When not in use, it is safe to store Touch on the charger.

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including but not limited to the following.

READ ALL INSTRUCTIONS BEFORE USE

DANGER—To reduce the risk of electrocution:

1. **Do not** place or store the Touch charger where it can fall or be pulled into a tub or sink.
2. **Do not** place the charger in or drop the charger into water or other liquid, or use the charger while bathing.
3. **Do not** reach for a charger that has fallen into water. Unplug immediately.

WARNING – to reduce the risk of burns, electrocution, fire or injury to persons

1. Close supervision is necessary when this product is used by, on, or near children or physically challenged persons to avoid improper application.
2. Use this product only for its intended use as described in this manual. Do not use attachments or products not recommended by the manufacturer.
3. Keep the cord away from heated surfaces.
4. Never use while sleepy or drowsy.
5. Do not use outdoors or where oxygen is being administered.

Never operate the Touch if it is damaged or if the applicator tip is missing or damaged. Do not operate the charger if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Call Clarisonic at 888-525-2747.

SAVE THESE INSTRUCTIONS

The Touch applicator, dispensing cap and charger have no serviceable parts. The supply cord cannot be replaced. If the cord is damaged, the charger must be unplugged and replaced. Call Clarisonic at 888-525-2747.

NOTE: This equipment has been tested and found to comply with Part 18 of the FCC Rules. Refer to our website for further information.

IMPORTANT CAUTIONS AND WARNINGS

If you are under a doctor's care for, or if you think you may have a skin condition or your eyes are sensitive to skin care products, please consult your physician before using the Clarisonic Touch. Refer to our website for more information.

We recommend that you **do not** share applicator tips if you and/or any other Clarisonic user suffer from any skin conditions. **Do not** use the Touch directly on the eye, the eyelid, or within the soft area around your eye. See [show illustration of okay/not okay area] [other areas to specifically avoid?]

The Touch and the Sea Serum are optimized to work synergistically for optimal results. We recommend that you **do not use** the Touch with other products. **Do not** use the Touch with any product containing mechanical exfoliants such as microbeads and shells. PBL is not responsible for any adverse results from using the Touch with any product other than the included Sea Serum.

See and Feel Guarantee

We are happy to offer a 90 day money back guarantee. If you are for any reason dissatisfied with the purchase of your Clarisonic Touch from an authorized dealer, call Customer Service at 888-525-2747 to make arrangements for the return.

In addition, every Clarisonic Touch comes with a two-year warranty covering any problems resulting from everyday use of the product. Call Customer Service at 888-525-2747 if you feel your Clarisonic is not operating properly.

For more comprehensive information, new products, product user guide updates, and to stay connected, visit **clarisonic.com**.

© Pacific Bioscience Laboratories, Inc. 2004-09.

All rights reserved. Clarisonic and [product name/s] are trademarks/registered trademarks of Pacific Bioscience Laboratories, Inc.

Patent listing:

Pacific Bioscience Laboratories, Inc.

13222 SE 30th Street #A1

Bellevue WA 98005

Toll Free: 888-5-CLARISONIC or 888-525-2747

Assembled in USA

Registration information

Register your Clarisonic Touch...

Also available from Clarisonic...

Touch complements sonic cleansing: The Clarisonic cleansing systems gently cleanse the skin and prepare it to better absorb serums and moisturizers. The skin care ritual continues with Touch's effective application of anti-aging serum.

Clarisonic Firming Sea Serum

[Note: the serum copy will morph into its own user documentation for use in the Touch and the replenishment packaging. Copy is not yet final]

Rich in antioxidants, Clarisonic Firming Sea Serum provides a multi-layered approach to balance and protect the skin. Nutrient and polysaccharide-rich marine and botanical extracts create a skin firming complex that, with sonic application, results in skin that is immediately smoother, tighter and refreshed.

Benefits include:

- Moisture protection for immediate and long-lasting hydration
- Promotes cellular metabolism and collagen synthesis (in vitro)
- Reduces the appearance of fine lines and wrinkles
- Soothes and calms irritation-prone skin
- Brightens and evens skin tone
- Targets dark under-eye circles

Key ingredients:

Canadian Willow Herb

Oat Kernel Extract

Hydrolyzed Elastin and Soluble Collagen

Kigelia Africana Fruit Extract

Glycosaminoglycans

Prickly Pear

Vitamins C and B3

Gentle for all skin types. Paraben and Phthalate Free. No animal testing. Developed and made in USA.

Sea Serum Ingredients

Water (Aqua), Epilobium Angustifolium Leaf/Stem Extract, Glycerin, Hydrolyzed Glycosaminoglycans, PEG-8 Dimethicone, Kigelia Africana Fruit Extract, Silanetriol Trehalose Saccharomyces Ferment, Avena Sativa (Oat) Kernel Extract, Ceratonia Siliqua Gum, Laminaria Saccharina Extract, Algae Extract, Pullulan, Skeletonema Costatum Amino Arginine, Arginine PCA, Plankton Extract, Soluble Collagen, Chondrus Extract, Chondrus Crispus

(Carrageenan), Niacinamide, Ascorbic Acid, Bisabolol, (Ginger) Root Extract, Opuntia Coccinellifera Fruit Extract, Sodium PCA, Urea, Hyaluronate, Polyquaternium-51, Ethyl Lauroyl Arginate HCL, Butylene Glycol, Acrylate Crosspolymer, Chlorphenesin, Disodium EDTA, Sodium Hydroxide