



# TRACE 1300 and TRACE 1310

Gas Chromatographs

## Hardware Manual

PN 31715002 Revision F January 2016

**Thermo**  
SCIENTIFIC

© 2016 Thermo Fisher Scientific Inc. All rights reserved.

TRACE 1300 and TRACE 1310 are trademarks of Thermo Fisher Scientific. Microsoft® is a registered trademark of Microsoft. Adobe® is a registered trademark of Adobe Systems Incorporated in the United States and/or other countries. Swagelok® is a registered trademark of Swagelok Company. Viton® and Vespel® are registered trademarks of DuPont. SilFlow®, NoVent®, and FingerTite® are registered trademarks of SGE Analytical Science. All other trademarks are the property of Thermo Fisher Scientific and its subsidiaries.

Published by Thermo Fisher Scientific S.p.A., Strada Rivoltana 20090 Rodano-Milan, Italy  
Tel: +39 02 95059373; Fax: +39 02 95059388

Thermo Fisher Scientific Inc. provides this document to its customers with a product purchase to use in the product operation. This document is copyright protected and any reproduction of the whole or any part of this document is strictly prohibited, except with the written authorization of Thermo Fisher Scientific Inc.

The contents of this document are subject to change without notice. All technical information in this document is for reference purposes only. System configurations and specifications in this document supersede all previous information received by the purchaser.

**Thermo Fisher Scientific Inc. makes no representations that this document is complete, accurate or error-free and assumes no responsibility and will not be liable for any errors, omissions, damage or loss that might result from any use of this document, even if the information in the document is followed properly.**

This document is not part of any sales contract between Thermo Fisher Scientific Inc. and a purchaser. This document shall in no way govern or modify any Terms and Conditions of Sale, which Terms and Conditions of Sale shall govern all conflicting information between the two documents.

Release history:

First Edition, released March 2012 - *“Original Instructions”*

Second Edition September 2012; Third Edition October 2013; Fourth Edition January 2014, Fifth Edition March 2015; Sixth Edition, January 2016

**For Research Use Only. Not for use in diagnostic procedures.**

## TRACE 1300 and TRACE 1310 Hardware Manual, PN 31715002, Revision F

|   | Strongly Agree | Agree | Neutral | Disagree | Strongly Disagree |
|---|----------------|-------|---------|----------|-------------------|
| The manual is well organized.                       | 1              | 2     | 3       | 4        | 5                 |
| The manual is clearly written.                      | 1              | 2     | 3       | 4        | 5                 |
| The manual contains all the information I need.     | 1              | 2     | 3       | 4        | 5                 |
| The instructions are easy to follow.                | 1              | 2     | 3       | 4        | 5                 |
| The instructions are complete.                      | 1              | 2     | 3       | 4        | 5                 |
| The technical information is easy to understand.    | 1              | 2     | 3       | 4        | 5                 |
| Examples of operation are clear and useful.         | 1              | 2     | 3       | 4        | 5                 |
| The figures are helpful.                            | 1              | 2     | 3       | 4        | 5                 |
| I was able to operate the system using this manual. | 1              | 2     | 3       | 4        | 5                 |

fold — — If not, please comment below. Attach additional sheets if necessary.

|  |  |
|--|--|
|  |  |
|  |  |
|  |  |
|  |  |
|  |  |
|  |  |
|  |  |

## Customer Registration Card

Register now...and receive all the privileges associated with being a Thermo Fisher Scientific product user including customer support, application reports, and technical reports.

### MY ORGANIZATION IS: (Check only one)

- ☐ Commercial (for profit) lab
- ☐ Government lab
- ☐ Hospital/Clinic
- ☐ Industrial lab
- ☐ Research Institute
- ☐ University/College
- ☐ Veterinary
- ☐ Other \_\_\_\_\_

### MY PRIMARY APPLICATION IS: (Check only one)

- ☐ Analytical
- ☐ Biomedical
- ☐ Clinical/Toxicology
- ☐ Energy
- ☐ Environmental
- ☐ Food/Agricultural
- ☐ Forensic/Toxicology
- ☐ Pharmaceutical
- ☐ Research/Education
- ☐ Other \_\_\_\_\_

### MY PRIMARY JOB FUNCTION IS: (Check only one)

- ☐ Administration
- ☐ Lab management
- ☐ Operator
- ☐ Other \_\_\_\_\_

fold — —

|                     |                      |
|---------------------|----------------------|
| Name _____          | Title _____          |
| Company _____       |                      |
| Address _____       |                      |
| City/State _____    | Postal Code _____    |
| Country _____       |                      |
| Telephone _____     | Ext. _____           |
| Serial Number _____ | Date purchased _____ |

### Fold and mail or e-mail to:

Editor, Technical Publications  
Thermo Fisher Scientific S.p.A.  
Strada Rivoltana km 4  
20090 Rodano (MI)  
Italy

Editor, Technical Publications  
Thermo Fisher Scientific GC-GC/MS  
2215 Grand Avenue Parkway  
Austin TX 78728-3812  
Unites States of America





## Declaration

Manufacturer: **Thermo Fisher Scientific**

Thermo Fisher Scientific is the manufacturer of the instrument described in this manual and, as such, is responsible for the instrument safety, reliability and performance only if:

- installation
- re-calibration
- changes and repairs

have been carried out by authorized personnel and if:

- the local installation complies with local law regulations
- the instrument is used according to the instructions provided and if its operation is only entrusted to qualified trained personnel

Thermo Fisher Scientific is not liable for any damages derived from the non-compliance with the aforementioned recommendations.

**Thermo Fisher Scientific S.p.A.**

Strada Rivoltana, 20090 Rodano - Milan - Italy — Tel: +39 02 950591 - Fax: +39 02 9505276

## Regulatory Compliance

Thermo Fisher Scientific performs complete testing and evaluation of its products to ensure full compliance with applicable domestic and international regulations.

When the system is delivered to you, it meets all pertinent electromagnetic compatibility (EMC) and safety standards.

### Safety

This device complies with the following safety standards according to Machinery Directive 2006/42/EC and Low Voltage Directive 2006/95/EC.

- International Electrotechnical Commission (IEC): 61010-1:2001 (2nd edition) - 61010-2-010:2003 (2nd edition) - 61010-2-081:2001 (1st edition) + A1:(2003)
- National differences: CAN/CSA C22.2 No. 61010-1 (2nd Edition) - UL 61010-1 (2nd Edition)
- EuroNorm (EN): 61010-1:2001 (2nd Edition) – 61010-2-010:2004 (2nd Edition) - 61010-2-081:2002 (1st Edition)

### Electromagnetic Compatibility

This device complies with the following regulations on Electromagnetic Compatibility (EMC) and Radio Frequency Interference (RFI) according to directive 2004/108/EC:

- CISPR 11/EN 55011: Group 1 Class A
- IEC/EN 61326-1:2012 (



**IMPORTANT:** Class A equipment is intended for use in an industrial environment. In others environments there may be potential difficulties in ensuring electromagnetic compatibility, due to the conducted as well as radiated disturbances.

## FCC Compliance Statement

THIS DEVICE COMPLIES WITH PART 15 OF THE FCC RULES. OPERATION IS SUBJECT TO THE FOLLOWING TWO CONDITIONS: (1) THIS DEVICE MAY NOT CAUSE HARMFUL INTERFERENCE, AND (2) THIS DEVICE MUST ACCEPT ANY INTERFERENCE RECEIVED, INCLUDING INTERFERENCE THAT MAY CAUSE UNDESIRE OPERATION.



**CAUTION** Read and understand the various precautionary notes, signs, and symbols contained inside this manual pertaining to the safe use and operation of this product before using the device.

## Notice on Lifting and Handling of Thermo Scientific Instruments

For your safety, and in compliance with international regulations, the physical handling of this Thermo Fisher Scientific instrument *requires a team effort* to lift and/or move the instrument. This instrument is too heavy and/or bulky for one person alone to handle safely.

## Notice on the Proper Use of Thermo Scientific Instruments

In compliance with international regulations: Use of this instrument in a manner not specified by Thermo Fisher Scientific could impair any protection provided by the instrument.

## Notice on the Susceptibility to Electromagnetic Transmissions

Do not use radio frequency transmitters, such as mobile phones, in close proximity to the instrument.

## WEEE Compliance

This product is required to comply with the European Union's Waste Electrical & Electronic Equipment (WEEE) Directive 2012/19/EU. It is marked with the following symbol:



Thermo Fisher Scientific has contracted with one or more recycling or disposal companies in each European Union (EU) Member State, and these companies should dispose of or recycle this product. See [www.thermoscientific.com/rohswcee](http://www.thermoscientific.com/rohswcee) for further information on Thermo Fisher Scientific's compliance with these Directives and the recyclers in your country.

## WEEE Konformität

Dieses Produkt muss die EU Waste Electrical & Electronic Equipment (WEEE) Richtlinie 2012/19/EU erfüllen. Das Produkt ist durch folgendes Symbol gekennzeichnet:



Thermo Fisher Scientific hat Vereinbarungen mit Verwertungs-/Entsorgungsfirmen in allen EU-Mitgliedsstaaten getroffen, damit dieses Produkt durch diese Firmen wiederverwertet oder entsorgt werden kann. Mehr Information über die Einhaltung dieser Anweisungen durch Thermo Fisher Scientific, über die Verwerter, und weitere Hinweise, die nützlich sind, um die Produkte zu identifizieren, die unter diese RoHS Anweisung fallen, finden sie unter [www.thermoscientific.com/rohswcee](http://www.thermoscientific.com/rohswcee).

## Conformité DEEE

Ce produit doit être conforme à la directive européenne (2012/19/EU) des Déchets d'Equipements Electriques et Electroniques (DEEE). Il est marqué par le symbole suivant:



Thermo Fisher Scientific s'est associé avec une ou plusieurs compagnies de recyclage dans chaque état membre de l'union européenne et ce produit devrait être collecté ou recyclé par celles-ci. Davantage d'informations sur la conformité de Thermo Fisher Scientific à ces directives, les recycleurs dans votre pays et les informations sur les produits Thermo Fisher Scientific qui peuvent aider la détection des substances sujettes à la directive RoHS sont disponibles sur [www.thermoscientific.com/rohswcee](http://www.thermoscientific.com/rohswcee).

## Conformità RAEE

Questo prodotto è marcato con il seguente simbolo in conformità alla direttiva europea 2012/19/EU (RAEE) sui rifiuti di apparecchiature elettriche ed elettroniche:



Thermo Fisher Scientific si è accordata con una o più società di riciclaggio in ciascun Stato Membro della Unione Europea (EU), e queste società dovranno smaltire o riciclare questo prodotto. Per maggiori informazioni vedere il sito [www.thermoscientific.com/rohsweee](http://www.thermoscientific.com/rohsweee).

## Conformidad RAEE

Este producto es marcado con el siguiente símbolo en conformidad a la Directiva 2012/19/EU de la Unión Europea sobre los residuos de aparatos eléctricos y electrónicos:



Thermo Fisher Scientific ha contratado una o más empresas de reciclaje para tratar residuos en cada Estado Miembro de la Unión Europea, y estas empresas deberían reciclar o eliminar este producto. Referirse a [www.thermoscientific.com/rohsweee](http://www.thermoscientific.com/rohsweee) para una mayor información sobre la conformidad de Thermo Fisher Scientific con estas Directivas y para las empresas de reciclaje en su país.

# Contents

|  |             |
|--|-------------|
| <b>Preface</b>   | <b>xvii</b> |
| About Your System                                      | xviii       |
| Power Rating   | xviii       |
| Contacting Us  | xix         |
| Related Documentation                                  | xix         |
| Safety Alerts and Important Information                | xix         |
| Special Notices  | xix         |
| Safety Symbols and Signal Words                        | xx          |
| Instrument Markings and Symbols                        | xxi         |
| Hydrogen Safety Precautions                            | xxii        |
| Using Hydrogen with TRACE 1300/TRACE 1310              | xxiii       |
| Hydrogen Connection Guidelines                         | xxiv        |
| Purchasing Hydrogen                                    | xxvi        |
| Properly Storing Hydrogen                              | xxvi        |
| Hydrogen Safety Codes, Standards and References        | xxvii       |
| Hazardous Substances Precautions                       | xxix        |
| Venting Toxic Gases                                    | xxix        |
| Liquid Nitrogen Safety Precautions                     | xxix        |
| Carbon Dioxide Safety Precautions                      | xxx         |
| <br><b>Chapter 1</b>                                   |             |
| <b>Installation</b>                                    | <b>1</b>    |
| Positioning the TRACE 1300/TRACE 1310                  | 3           |
| Installing the External Accessories                    | 3           |
| Installing the Injector and Detector Modules           | 4           |
| Installing an Injector Module                          | 4           |
| Installing a Detector Module                           | 6           |
| Making the Gas Supply Plumbing Connections             | 7           |
| Testing for Leaks                                      | 12          |
| Connecting the Oven Cryogenic System                   | 13          |
| Oven Cryo System with Carbon Dioxide Connection        | 13          |
| Oven Cryo System with Liquid Nitrogen Connection       | 16          |
| Connecting the PTV/PTVBKF Cryogenic System             | 18          |
| PTV/PTVBKF Cryo System with Carbon Dioxide Connection  | 18          |
| PTV/PTVBKF Cryo System with Liquid Nitrogen Connection | 21          |

|   |               |
|---|---------------|
| Coupling to a Mass Spectrometer . . . . .   | 23            |
| Making the Duct to Couple with an ISQ Series, TSQ 8000 Series,<br>DSQ II, or ITQ Mass Spectrometer. . . . . | 23            |
| Making the Duct to Couple with a TSQ Quantum Mass<br>Spectrometer . . . . .                                 | 26            |
| Making the Duct to Couple with a DFS, IRMS, or ICP-MS Mass<br>Spectrometer . . . . .                        | 29            |
| Setting Handshake Parameters. . . . .   | 31            |
| Installing the Autosampler. . . . .   | 32            |
| Mounting an Autosampler on the GC. . . . .  | 32            |
| Connecting the Autosampler. . . . .   | 34            |
| Connect the Autosampler to the GC System. . . . .   | 37            |
| Setting Autosampler Handshake Parameters . . . . .  | 41            |
| Installing the Data System Software . . . . .   | 42            |
| Making Power Connections . . . . .  | 43            |
| Setting the LAN Communication . . . . .   | 45            |
| Making the LAN Setup. . . . .   | 46            |
| Configuring the Data System . . . . .   | 50            |
| Column Installation Requirements . . . . .  | 50            |
| Using the Correct Fittings. . . . .   | 50            |
| Installing the Adapters for Encapsulated Graphite Ferrules . . . . .  | 52            |
| Installing the Column Rack. . . . .   | 56            |
| Installing the Column the First Time . . . . .  | 57            |
| <br><b>Chapter 2   Performing Routine Maintenance . . . . .</b>   | <br><b>65</b> |
| Read Me First . . . . .   | 66            |
| Maintenance Supplies and Tools. . . . .   | 68            |
| Cleaning Stainless Steel Components . . . . .   | 68            |
| Maintenance Button . . . . .  | 69            |
| Powering On the TRACE 1300/TRACE 1310. . . . .  | 69            |
| Shutting Down the TRACE 1300/TRACE 1310 . . . . .   | 71            |
| Cleaning the Instrument Externally. . . . .   | 72            |
| Replacing a Column . . . . .  | 73            |
| <br><b>Chapter 3   Performing Injectors Routine Maintenance. . . . .</b>                                    | <br><b>79</b> |
| Maintaining a Split/Splitless Injector (SSL). . . . .   | 80            |
| Replacing the SSL Septum. . . . .   | 83            |
| Cleaning or Replacing the SSL Glass Liner . . . . .   | 84            |
| Replacing a SSL Broken Liner . . . . .  | 86            |
| Replacing the SSL Carrier and Split Lines Filters. . . . .  | 88            |
| Replacing the SSL Body Head O-Rings. . . . .  | 89            |

|   |            |
|---|------------|
| Maintaining a Split/Splitless Injector with Backflush (SSLBKF) . . . . .                        | 91         |
| Replacing the SSLBKF Septum . . . . .   | 94         |
| Cleaning or Replacing the SSLBKF Glass Liner. . . . .   | 95         |
| Replacing a SSLBKF Broken Liner . . . . .   | 97         |
| Replacing the SSLBKF Carrier and Split Lines Filters . . . . .                                  | 99         |
| Replacing the SSLBKF Body Head O-Rings . . . . .  | 100        |
| Maintaining a Gas Sampling Valve Injector (GSV) . . . . .                                       | 102        |
| Connecting the Sample In and Out Lines . . . . .  | 103        |
| Replacing the Carrier and Split Lines Filters . . . . .   | 105        |
| . . . . .   | 105        |
| Replacing the Sample Loop . . . . .   | 106        |
| Maintaining an Instant Connect Helium Saver Injector (HeS-S/SL) . . . . .                       | 108        |
| Replacing the Septum . . . . .  | 111        |
| Cleaning or Replacing the Glass Liner. . . . .  | 112        |
| Replacing the Carrier and Split Lines Filters . . . . .   | 114        |
| Replacing the Body Head O-Rings . . . . .   | 115        |
| Maintaining a Programmable Temperature Vaporizing Injector (PTV) . . . . .                      | 117        |
| Replacing the PTV Septum . . . . .  | 120        |
| Cleaning or Replacing the PTV Glass Liner . . . . .   | 121        |
| Replacing the PTV Broken Liner. . . . .   | 123        |
| Replacing the PTV Carrier and Split Lines Filters . . . . .                                     | 125        |
| Maintaining a Programmable Temperature Vaporizing Injector with<br>Backflush (PTVBKF) . . . . . | 127        |
| Replacing the PTVBKF Septum . . . . .   | 130        |
| Cleaning or Replacing the PTVBKF Glass Liner . . . . .  | 131        |
| Replacing the PTVBKF Broken Liner . . . . .   | 133        |
| Replacing the PTVBKF Carrier and Split Lines Filters . . . . .                                  | 135        |
| <b>Chapter 4 Performing Detectors Routine Maintenance . . . . .</b>                             | <b>137</b> |
| Maintaining a Flame Ionization Detector (FID) . . . . .   | 138        |
| Cleaning or Replacing the FID Collecting Electrode . . . . .                                    | 141        |
| Replacing the FID Ignition Glow-plug . . . . .  | 144        |
| Maintaining a Nitrogen Phosphorous Detector (NPD) . . . . .                                     | 147        |
| Replacing the NPD Thermionic Source . . . . .   | 150        |
| Cleaning or Replacing the NPD Collecting Electrode . . . . .                                    | 156        |
| Maintaining a Thermal Conductivity Detector (TCD) . . . . .                                     | 163        |
| Bake-out Procedure . . . . .  | 164        |
| Measuring the Carrier Gas Flow Rate . . . . .   | 165        |
| Shutting Down the TCD . . . . .   | 165        |
| Maintaining an Electron Capture Detector (ECD) . . . . .  | 166        |
| Detector Chemical Contamination . . . . .   | 168        |
| Wipe Test . . . . .   | 169        |
| Cleaning or Replacing the Collecting Electrode (Anode) . . . . .                                | 169        |



|                  |  |            |
|------------------|--|------------|
|                  | Maintaining a Flame Photometric Detector (FPD) .....                                   | 172        |
|                  | Installing the FPD Detector .....  | 176        |
|                  | Removing the FPD Detector .....  | 178        |
|                  | Cleaning or Replacing the FPD Jet .....  | 179        |
|                  | Cleaning or Replacing the FPD Interferential Filter .....                              | 183        |
|                  | Replacing the FPD Ignition Glow-plug .....   | 187        |
|                  | Maintaining a Pulsed Discharge Detector (PDD) .....                                    | 193        |
| <b>Chapter 5</b> | <b>GC Main Frame Advanced Maintenance .....</b>  | <b>195</b> |
|                  | Removing/Replacing the GC Top Cover .....  | 196        |
|                  | Removing/Replacing the GC Left Side Panel .....  | 197        |
|                  | Removing/Replacing the GC Right Side Panel .....                                       | 199        |
|                  | Removing/Replacing the GC Back Cover .....   | 201        |
|                  | Removing/Replacing the GC Front Door Cover .....                                       | 203        |
|                  | Removing/Replacing the Electronic Module .....   | 208        |
|                  | Replacing the Oven Heater Baffle .....   | 211        |
|                  | Replacing the Oven Heater Temperature Sensor .....                                     | 215        |
|                  | Replacing the Oven Motor .....   | 219        |
|                  | Replacing the Flap Motor .....   | 225        |
| <b>Chapter 6</b> | <b>Injectors Advanced Maintenance .....</b>  | <b>229</b> |
|                  | Baking-out Contaminants from SSL, SSLBKF, HeS-S/SL, PTV, and<br>PTVBKF Injectors. .... | 230        |
|                  | Removing/Replacing an Injector Module .....  | 231        |
|                  | Cleaning the SSL Injector Body .....   | 234        |
|                  | Cleaning the SSLBKF Injector Body .....  | 238        |
|                  | Cleaning the HeS-S/SL Injector Body .....  | 242        |
|                  | Cleaning the PTV Injector Head Assembly .....  | 246        |
|                  | Cleaning the PTVBKF Injector Head Assembly .....                                       | 253        |
| <b>Chapter 7</b> | <b>Detectors Advanced Maintenance .....</b>  | <b>261</b> |
|                  | Removing/Replacing a Detector Module .....   | 262        |
|                  | Measuring the FID Gas Flows .....  | 265        |
|                  | Cleaning or Replacing the FID Jet .....  | 271        |
|                  | Measuring the NPD Gas Flows .....  | 276        |
|                  | Cleaning or Replacing the NPD Jet .....  | 285        |
|                  | Measuring the FPD Gas Flows .....  | 293        |
|                  | Cleaning or Replacing the FPD Mirror Metal Plug .....                                  | 295        |
|                  | Cleaning or Replacing the FPD Filter-side Heat Shields .....                           | 300        |
|                  | Cleaning or Replacing the FPD Flame-side Heat Shields .....                            | 308        |
|                  | Replacing the FPD Photomultiplier Tube .....   | 317        |
| <b>Chapter 8</b> | <b>Installing Optional Kits .....</b>  | <b>325</b> |
|                  | Installing the Oven Exhaust Kit .....  | 326        |

|  |                |
|--|----------------|
| Installing the Merlin Microseal High Pressure Valve Kit . . . . .              | 327            |
| Introduction . . . . .   | 327            |
| Getting Started . . . . .  | 327            |
| Installing the Purge & Trap Adapter Kit on the SSL/SSLBKF Injector . . . . .   | 332            |
| Installing the Packed Column Adapters . . . . .                                | 338            |
| Introduction . . . . .   | 338            |
| Getting Started . . . . .  | 338            |
| Installing the HS Adapter Kit on the SSL/SSLBKF Injector . . . . .             | 344            |
| Installing the Large Volume Splitless Kit . . . . .                            | 354            |
| Installing the Manual On/Off Valve for Single Gas Line . . . . .               | 357            |
| Connecting a SSL/PTV Backflush System . . . . .                                | 360            |
| Connecting a SSL/PTV Backflush System for High Temperature . . . . .           | 367            |
| Connecting a GSV Backflush System . . . . .                                    | 373            |
| Installing the NoVent Microfluidics . . . . .                                  | 376            |
| Connecting the NoVent Microfluidics Module to the TRACE<br>1300/1310 . . . . . | 376            |
| Installing the Mounting Bracket . . . . .                                      | 379            |
| Preparing the NoVent Microfluidics Restrictor Tubing . . . . .                 | 381            |
| Attaching the Ferrule and Nut to the GC Column . . . . .                       | 384            |
| Attaching the New Tubing to the Transfer Line . . . . .                        | 385            |
| Connecting the Capillaries to the Microfluidics Splitter . . . . .             | 390            |
| Configuring the Post-Column . . . . .  | 391            |
| Using the Module . . . . .   | 394            |
| Installing a FTIR Make-up Module . . . . .                                     | 396            |
| Installing the Hot Injection Adapter Kit on the SSL/SSLBKF Injector . . . . .  | 400            |
| Performing the Dual FPD Detector Configuration . . . . .                       | 409            |
| <br><b>Chapter 9    Adding Modules . . . . .</b>                               | <br><b>421</b> |
| Adding a SSL, SSLBKF, PTV, or PTVBKF Injector Module . . . . .                 | 422            |
| Adding a GSV Injector Module . . . . .   | 425            |
| Adding a FID, TCD/TCD In-Series, ECD, or FPD Detector Module . . . . .         | 431            |
| Adding a NPD Detector Module . . . . .   | 435            |
| Adding an Aux Temperature/Cryo Module . . . . .                                | 443            |
| Adding a Helium Saver Injector Module . . . . .                                | 448            |
| Installing the Module . . . . .  | 448            |
| Installing the Column . . . . .  | 454            |
| Checking for Leaks . . . . .   | 456            |

|   |     |
|---|-----|
| Adding a PDD Module . . . . .   | 458 |
| Getting Started . . . . .   | 458 |
| Installing a PDD Module . . . . .                                       | 459 |
| Plumbing the Gas Lines. . . . .   | 463 |
| Connecting the High Voltage and Pulses Cables . . . . .                 | 465 |
| Installing the Capillary Column . . . . .                               | 470 |
| Installing the Packed Column . . . . .                                  | 470 |
| Testing for Leaks . . . . .   | 470 |
| Performing Initial Power Up . . . . .                                   | 471 |
| Adding a Generic Detector Interface . . . . .                           | 472 |
| Preliminary Operations . . . . .  | 472 |
| Getting Started . . . . .   | 472 |
| Removing the GC Back Cover . . . . .                                    | 473 |
| Assembling the GDI Electrical Interface . . . . .                       | 474 |
| Replacing the Encapsulated Flow Restrictors . . . . .                   | 476 |
| Installing and Connecting the GDI Electrical Interface . . . . .        | 477 |
| Connecting a GDI Interface to the TRACE 1310 Auxiliary Oven. . . . .    | 481 |
| Installing a GDI Mechanical Module . . . . .                            | 482 |
| Connecting the Detector Gas Tubings to the Manifolds . . . . .          | 484 |
| Connecting Heater and Signal Cables . . . . .                           | 488 |
| Restarting the GC . . . . .   | 489 |
| Performing the Third-party Detector Start-up and Optimization . . . . . | 489 |
| Configuring and Setting GDI Detector . . . . .                          | 490 |
| Adding an Analog Output Interface . . . . .                             | 491 |

## **Chapter 10 Adding Systems . . . . .495**

|  |     |
|--|-----|
| Adding the Oven Cryo System . . . . .  | 496 |
| Oven Cryo System Overview. . . . .   | 497 |
| Installing the Oven Cryo System . . . . .                                      | 499 |
| Adding the PTV and PTVBKF Cryo System. . . . .                                 | 510 |
| PTV/PTVBKF Cryo System Overview . . . . .                                      | 511 |
| Installing the PTV/PTVBKF Cryo System . . . . .                                | 513 |
| Adding an Auxiliary Gas System . . . . .                                       | 521 |
| Auxiliary Gas Module Overview . . . . .  | 522 |
| Auxiliary Gas Interface Overview. . . . .                                      | 522 |
| Preliminary Operations . . . . .   | 524 |
| Getting Started . . . . .  | 525 |
| Installing the Auxiliary Gas Interface on the Left Wall of the Oven. . . . .   | 525 |
| Installing the Auxiliary Gas Interface on the Right Wall of the Oven . . . . . | 528 |
| Installing the Auxiliary Gas Interface on the Oven for HRMS. . . . .           | 532 |
| Installing and Connecting the Auxiliary Gas Module . . . . .                   | 534 |
| Adding the Hydrogen Sensor . . . . .   | 542 |

## **Chapter 11 Upgrade Equipment. . . . .549**

|  |     |
|--|-----|
| Upgrading a TRACE 1300 to a TRACE 1310 . . . . . | 550 |
|--|-----|

|  |            |
|--|------------|
| Upgrading a Stand Alone TRACE 1300/TRACE 1310 to MS Version. . . . .       | 556        |
| Preliminary Operations . . . . .   | 556        |
| Getting Started . . . . .  | 556        |
| Coupling with the ISQ Series / TSQ 8000 Series Mass Spectrometer . . . . . | 557        |
| Coupling with the DSQ II Mass Spectrometer . . . . .                       | 561        |
| Coupling with the ITQ Mass Spectrometer. . . . .                           | 566        |
| Coupling with the TSQ Quantum Mass Spectrometer . . . . .                  | 571        |
| Updating HMI Software From USB Stick. . . . .                              | 577        |
| <b>Chapter 12 Troubleshooting. . . . .</b>                                 | <b>581</b> |
| Investigating Power Supply Issues . . . . .                                | 582        |
| TRACE 1300/TRACE 1310 will not power-on . . . . .                          | 582        |
| Investigating Communication Issues . . . . .                               | 582        |
| Software is not communicating with the TRACE 1300/TRACE 1310 . . . . .     | 582        |
| TRACE 1300/TRACE 1310 does not start or is not ready . . . . .             | 582        |
| Cannot download methods to the TRACE 1300/TRACE 1310 . . . . .             | 582        |
| Sample data are not acquired. . . . .                                      | 583        |
| GC is not communicating with the PC. . . . .                               | 583        |
| Autosampler is not communicating with the PC . . . . .                     | 583        |
| Investigating Sensitivity Issues . . . . .                                 | 583        |
| Poor sensitivity or sudden loss in sensitivity. . . . .                    | 583        |
| Error Messages . . . . .   | 584        |
| TRACE 1300 Error Messages . . . . .  | 584        |
| TRACE 1310 Error Messages . . . . .  | 584        |
| Contacting Technical Support. . . . .                                      | 589        |
| <b>Glossary . . . . .</b>  | <b>591</b> |



# Preface

This manual contains detailed information about installing, maintaining, and troubleshooting the TRACE 1300/TRACE 1310 gas chromatograph.

This manual is organized as follows:

**Chapter 1, “Installation,”** provides instructions for installing and connecting the TRACE 1300/TRACE 1310 system.

**Chapter 2, “Performing Routine Maintenance,”** provides instructions for performing routine maintenance on the TRACE 1300 and TRACE 1310 modules.

**Chapter 3, “Performing Injectors Routine Maintenance,”** provides instructions for performing routine maintenance on the TRACE 1300 and TRACE 1310 injector modules.

**Chapter 4, “Performing Detectors Routine Maintenance,”** provides instructions for performing routine maintenance on the TRACE 1300 and TRACE 1310 detector modules.

**Chapter 5, “GC Main Frame Advanced Maintenance,”** describes TRACE 1300/TRACE 1310 components that do not require routine maintenance, but they need to be removed or replaced.

**Chapter 6, “Injectors Advanced Maintenance,”** describes TRACE 1300/TRACE 1310 injector modules that do not require routine maintenance, but troubleshooting may indicate they need to be cleaned or replaced.

**Chapter 7, “Detectors Advanced Maintenance,”** describes TRACE 1300/TRACE 1310 detector modules that do not require routine maintenance, but troubleshooting may indicate they need to be cleaned or replaced.

**Chapter 8, “Installing Optional Kits,”** describes how to install the optional kits available for the TRACE 1300/TRACE 1310.

**Chapter 9, “Adding Modules,”** describes how to install any added module that is available for the TRACE 1300/TRACE 1310.

**Chapter 10, “Adding Systems,”** describes how to install any added system that is available for the TRACE 1300/TRACE 1310. See the *TRACE 1300 and TRACE 1310 Spare Parts Guide* for information about ordering the equipment in this chapter.

Chapter 11, “Upgrade Equipment,” describes how to install any upgraded equipment that is available for the TRACE 1300/TRACE 1310.

Chapter 12, “Troubleshooting,” describes the symptom, and the remedy for each known issue with the TRACE 1300/TRACE 1310 gas chromatograph.

“Glossary,” contains definitions of terms used in this guide. It also includes abbreviations, acronyms, metric prefixes, and symbols.

## About Your System

Thermo Scientific systems provide high-caliber gas chromatography (GC) instrumentation. Your TRACE 1300/TRACE 1310 GC system can be a stand-alone unit or coupled with other instruments. GC represents a powerful analytical separation technique. Complex mixtures of individual compounds can be injected into the GC, either manually or by using an autosampler, and then separated the eluate for presentation to the detector. The detector generates signals of the GC eluate and its components. These signals are then processed by a Thermo Scientific Chromatography Data System for qualitative identification, as well as accurate and precise quantification of the individual compounds present in the sample.

**IMPORTANT** Thermo Scientific systems optimize the separation and detection capabilities of GC by providing high performance analytical capabilities for both research, and routine applications. More information about the use of this system can be found in related documentation sources, and by using the provided contact information.



**WARNING** Thermo Scientific systems operate safely and reliably under carefully controlled environmental conditions. If the equipment is used in a manner not specified by the manufacturer, the protections provided by the equipment might be impaired. If you maintain a system outside the specifications listed in this guide, failures of many types, including personal injury or death, might occur. The repair of instrument failures caused by operation in a manner not specified by the manufacturer is specifically excluded from the standard warranty and service contract coverage.



**WARNING** Operation of this system requires the use of chemical substances with different hazard specifications. Before using any chemicals, read the hazard indications and information reported in the Safety Sheet supplied by the manufacturer, referring to the relevant CAS (Chemical Abstract Service) number.

## Power Rating

TRACE 1300/TRACE 1310 gas chromatograph

- 120 Vac  $\pm$ 10%, 50/60 Hz, 2000 VA
- 230 Vac  $\pm$ 10%, 50/60 Hz, 2000 VA

Detailed instrument specifications are in the Product Specifications or Product Brochure.



## Contacting Us

Thermo Fisher Scientific provides comprehensive technical assistance worldwide and is dedicated to the quality of our customer relationships and services.

Use <http://www.thermoscientific.com> address for products information.

Use <http://www.gc-gcms-customersupport.com/WebPage/Share/Default.aspx> address to contact your local Thermo Fisher Scientific office or affiliate GC-GC/MS Customer Support.

## Related Documentation

In addition to this guide, Thermo Scientific provides the following documents for the TRACE 1300 and TRACE 1310. These documents are also available on a “Print-By-Request” basis.

*TRACE 1300 and TRACE 1310 Document Set*, PN 31715000

- *TRACE 1300 and TRACE 1310 Preinstallation Requirements Guide*, PN 31715001
- *TRACE 1300 and TRACE 1310 Hardware Manual*, PN 31715002
- *TRACE 1300 and TRACE 1310 User Guide*, PN 31715003
- *TRACE 1300 and TRACE 1310 Spare Parts Guide*, PN 31715004

## Safety Alerts and Important Information

Make sure you follow the precautionary notices presented in this manual. The safety and other special notices appear in boxes.

## Special Notices

Notices includes the following:













**IMPORTANT** Highlights information necessary to prevent damage to software, loss of data, or invalid test results; or might contain information that is critical for optimal performance of the system.








**Note** Emphasizes important information about a task.

**Tip** Helpful information that can make a task easier.

## Safety Symbols and Signal Words

All safety symbols are followed by **WARNING** or **CAUTION**, which indicates the degree of risk for personal injury, instrument damage, or both. Cautions and warnings are following by a descriptor, such as **BURN HAZARD**. A **WARNING** is intended to prevent improper actions that could cause personal injury. Whereas, a **CAUTION** is intended to prevent improper actions that might cause personal injury, instrument damage, or both. You can find the following safety symbols on your instrument, or in this manual:

| Symbol  | Descriptor   |
|---|--|
|    | <b>BIOHAZARD:</b> Indicates that a biohazard <i>will, could, or might</i> occur.   |
|    | <b>BURN HAZARD:</b> Alerts you to the presence of a hot surface that <i>could</i> or <i>might</i> cause burn injuries.   |
|    | <b>ELECTRICAL SHOCK HAZARD:</b> Indicates that an electrical shock <i>could</i> or <i>might</i> occur.   |
|    | <b>FIRE HAZARD:</b> Indicates a risk of fire or flammability <i>could</i> or <i>might</i> occur.   |
|   | <b>EXPLOSION HAZARD.</b> Indicates an explosion hazard. This symbol indicates this risk <i>could</i> or <i>might</i> cause physical injury.  |
|  | <b>FLAMMABLE GAS HAZARD.</b> Alerts you to gases that are compressed, liquefied or dissolved under pressure and can ignite on contact with an ignition source. This symbol indicates this risk <i>could</i> or <i>might</i> cause physical injury. |
|  | <b>GLOVES REQUIRED:</b> Indicates that you must wear gloves when performing a task or physical injury <i>could</i> or <i>might</i> occur.  |
|  | <b>CLOTHING REQUIRED.</b> Indicates that you should wear a work clothing when performing a task or else physical injury <i>could</i> or <i>might</i> occur.  |
|  | <b>BOOTS REQUIRED.</b> Indicates that you must wear boots when performing a task or else physical injury <i>could</i> or <i>might</i> occur.   |
|  | <b>MATERIAL AND EYE HAZARD.</b> Indicates you must wear eye protection when performing a task.   |
|  | <b>HAND AND CHEMICAL HAZARD:</b> Indicates that chemical damage or physical injury <i>could</i> or <i>might</i> occur.   |
|  | <b>HARMFUL.</b> Indicates that the presence of harmful material <i>will, could, or might</i> occur.  |
|  | <b>INSTRUMENT DAMAGE:</b> Indicates that damage to the instrument or component <i>might</i> occur. This damage might not be covered under the standard warranty.   |

|   |   |
|---|---|
|  | <b>LIFTING HAZARD.</b> Indicates that a physical injury <i>could</i> or <i>might</i> occur if two or more people do not lift an object.   |
|  | <b>MATERIAL AND EYE HAZARD:</b> Indicates that eye damage <i>could</i> or <i>might</i> occur.   |
|  | <b>READ MANUAL:</b> Alerts you to carefully read your instrument's documentation to ensure your safety and the instrument's operational ability. Failing to carefully read the documentation <i>could</i> or <i>might</i> put you at risk for a physical injury.  |
|  | <b>TOXIC SUBSTANCES HAZARD:</b> Indicates that exposure to a toxic substance could occur and that exposure <i>could</i> or <i>might</i> cause personal injury or death.   |
|  | <b>LASER HAZARD.</b> Indicates that exposure to a laser beam <i>will</i> , <i>could</i> , or <i>might</i> cause personal injury.  |
|  | <b>RADIOACTIVE HAZARD.</b> Indicates that the presence of radioactive material <i>could</i> or <i>might</i> occur.  |
|  | For the prevention of personal injury, this general warning symbol precedes the <b>WARNING</b> safety alert word and meets the ISO 3864-2 standard. In the vocabulary of ANSI Z535 signs, this symbol indicates a possible personal injury hazard exists if the instrument is improperly used or if unsafe actions occur. This symbol and another appropriate safety symbol alerts you to an imminent or potential hazard that <i>could cause personal injury</i> . |













## Instrument Markings and Symbols

Table 1 explains the symbols used on Thermo Scientific instruments. Only a few of them are used on the TRACE 1300/1310. See the asterisk.

**Table 1.** Instrument Marking and Symbols (Sheet 1 of 2)

| Symbol | Description                         |
|--------|-------------------------------------|
| ===    | Direct Current                      |
| * ~    | Alternating Current                 |
| ≍      | Both direct and alternating current |
| 3 ~    | Three-phase alternating current     |
| ⏏      | Earth (ground) terminal             |
| ⏏      | Protective conductor terminal       |
| ⏏      | Frame or chassis terminal           |

**Table 1.** Instrument Marking and Symbols (Sheet 2 of 2)

| Symbol  | Description   |
|---|---|
|      | Equipotentiality  |
| *    | On (Supply)   |
| *    | Off (Supply)  |
|      | Equipment protected throughout by DOUBLE INSULATION or REINFORCED INSULATION (Equivalent to Class II of IEC 536)  |
| *    | Instruction manual symbol affixed to product. Indicates that the user must refer to the manual for specific <b>WARNING</b> or <b>CAUTION</b> information to avoid personal injury or damage to the product. |
|      | Caution, risk of electric shock   |
| *    | Caution, hot surface  |
| *    | Caution, biohazard  |
|    | In-position of a bistable push control  |
|    | Out-position of a bistable push control   |
| *  | Jack socket   |
| *  | Symbol in compliance to the Directive 2012/19/EU on Waste Electrical and Electronic Equipment (WEEE) placed on the European market after August, 13, 2005.  |

## Hydrogen Safety Precautions

Hydrogen is a colorless, odorless, highly flammable gas with the molecular formula  $H_2$  and an atomic weight of 1.00794, making it the lightest element. Hydrogen gas presents a hazard, as it is combustible over a wide range of concentrations; at ambient temperature and pressure, the range is from about 4 to 74.2% by volume.

Hydrogen has a flash point of  $-423\text{ }^{\circ}\text{F}$  ( $-253\text{ }^{\circ}\text{C}$ ) and an auto-ignition temperature of  $1040\text{ }^{\circ}\text{F}$  ( $560\text{ }^{\circ}\text{C}$ ). It has a very low ignition energy and the highest burning velocity of any gas. If hydrogen is allowed to expand rapidly from high pressure, it can self-ignite. Hydrogen burns with a flame that can be invisible in bright light.



**WARNING - EXPLOSION HAZARD** The use of hydrogen as a carrier gas is dangerous. Hydrogen is potentially explosive and must be used with extreme care. Any use of hydrogen gas must be reviewed by appropriate health and safety staff, and all installations of hydrogen systems must be performed to applicable codes and standards. Thermo Fisher Scientific assumes no liability for the improper use of hydrogen as a carrier gas.

Before you begin using hydrogen, conduct a risk assessment based on the quantity of hydrogen to be used and the conditions of your laboratory. Ask yourself:

“What hydrogen hazards associated with this project are most likely to occur?”

“What hydrogen hazards associated with this project have the potential to result in the worst consequences?”

- Try to reduce or eliminate the higher risks by using the proper ventilation to remove hydrogen gas before an ignitable concentration can accumulate. Also consider purging the hydrogen to further reduce hazards and ensure that anyone working with hydrogen has basic hydrogen safety training.
- As with laboratory safety in general, be sure to wear safety glasses, laboratory coats, gloves, and so on. Typically there are no specific requirements for gaseous hydrogen, other than eye protection when working with a compressed gas. If working with liquid (cryogenic) hydrogen, wear insulated gloves and protective shoes in addition to eye protection.
- Post “No Smoking” and “No Open Flames” signs to identify hydrogen sources and cylinders. Maintain, inspect, and leak-test all hydrogen sources regularly.
- Clearly mark all hydrogen shutoff valves and label permanent hydrogen piping as such at the supply or discharge point, and at regular intervals along its length. Where hydrogen gas piping passes through a wall, be sure to label both sides of the wall.
- Have contingency plans in place should an incident occur.
- Ensure that site emergency response team, as well as the local fire department, knows the location of all hydrogen storage tanks.

## Using Hydrogen with TRACE 1300/TRACE 1310

The use of hydrogen as a carrier gas, or as fuel gas for certain flame detectors, requires strict attention and compliance with special precautions due to the hazards involved.

**WARNING - EXPLOSION HAZARD** Hydrogen is a dangerous gas that, when mixed with air, could create an explosive mixture. The use of hydrogen as a carrier gas requires extreme caution. Special precautions must be taken because of the risk of explosion. When hydrogen is used as carrier gas the gas chromatograph must be equipped with a hydrogen sensor.



Never use hydrogen as carrier gas in your TRACE 1300/TRACE 1310 system unless your oven has a hydrogen sensor installed. Thermo Fisher Scientific FSEs are not authorized to install or repair any instrument using hydrogen as a carrier gas unless the instrument is equipped with the appropriate sensor.

If your oven does not have a hydrogen sensor already installed, contact your Thermo Fisher Scientific sales representative. To comply with instrument safety requirements, a Thermo Fisher Scientific FSE authorized service personnel should install the sensor into your TRACE 1300/TRACE 1310.

Hydrogen is a dangerous gas, particularly in an enclosed area when it reaches a concentration corresponding to its lower explosion level (4% in volume). An explosion hazard could develop in the oven when hydrogen is used as a carrier gas in the case oven elements are not perfectly connected to each other, or when the connection materials are worn out, broken, or otherwise faulty.

Use the following safety precautions when using hydrogen:

- Ensure that all hydrogen cylinders comply with the safety requirements for proper use and storage. Hydrogen cylinders and delivery systems must comply with local regulations.
- Make sure the gas supply is turned completely off when connecting hydrogen lines.
- Perform a leak test to ensure that the hydrogen lines are leak-tight before using the instrument. Repeat this test to eliminate all leaks.
- Ensure your TRACE 1300/TRACE 1310 has a Thermo Scientific hydrogen sensor installed for continuously monitoring the hydrogen level in the oven.

## Hydrogen Connection Guidelines

Use the following guidelines to safely connect hydrogen to your system:

- **Piping**—Hydrogen must be delivered to equipment using appropriate piping and be done in such a way as to pose essentially no hazard to end-users. Piping systems for the delivery of hydrogen should be designed and installed by a person qualified by specific training and experience with hydrogen piping systems.

Stainless steel is usually recommended because it is a safe, cost-effective material. Piping of *black iron* or copper must not be used, as the pipe can become brittle with age.

Elastomeric/plastic tubing of various plastics and polymers should not be used, unless the

tubing is approved for use with hydrogen. If elastomeric/plastic tubing is used for hydrogen gas delivery, the tubing should be tested for hydrogen permeability to minimize leakage.

The hydrogen piping system must be flexible enough to endure routine thermal expansion and contraction. The system should also include considerations for the most severe condition of temperature and pressure expected during service. Piping and supports must be able to withstand static loading introduced by such things as ice and snow; and dynamic loading from high wind and earthquake.

Caution should be used if burying hydrogen piping. Proper controls should be used to protect against damage and corrosion, and also to prevent Hydrogen from entering a building if there is any leakage.

- **Fittings**—All fittings must be of the proper type approved or designed for use with hydrogen gas. Use as few fittings as possible to minimize the potential for leaks. After installation, ensure that leak testing is carried out prior to system use, and on a regular basis.

There must be no PTFE tape or other things like *plumber's putty* used to enhance a seal, as this actually is a detriment to a good seal. Ideally the best installation would use stainless steel tubing with appropriate gas-tight fittings.

Welding is usually preferred for joints in hydrogen piping systems since welding provides a better connection and reduces the potential for leaks compared to mechanical fittings. Soft solder joints are not permitted for hydrogen systems (due to the low melting point of soft solder and its potential for brittle failure at cryogenic temperatures). Braze joints are permitted, but such joints should be protected against the possibility of external fire.

Tubing connections should be clamped to barbed or press-fit type connections. Hose clamps or *jubilee clamps* must not be used.

- **Valves**—All valves must be suitable for hydrogen service and for the specific operating conditions. Valves, including regulators, must not be used for hydrogen, unless they are designed and identified for such a use. Ball valves are often chosen because of their superior leak tightness through the valve seat. Pneumatic operators are usually chosen for remotely operated valves so that potential ignition sources (electricity) are remote from the valve.

Manual shutoff valves should be provided near each point of use, within immediate reach. If a hydrogen cylinder or hydrogen generation system is located within immediate reach, a separate point-of-use shutoff valve is usually not necessary.

Line regulators that have their source away from the point of use should have a manual shutoff valve near the point of use.

An emergency gas shutoff device in an accessible location outside the use area should be provided in addition to the manual point-of-use valve in each educational and instructional laboratory space that has a piped gas supply system.



If necessary, the piping system should have uninterruptible pressure relief. The pressure relief system should be designed to provide a discharge rate sufficient to avoid further pressure increase and should vent to a safe location outside or to a ventilation system exhaust.

## Purchasing Hydrogen

Use the following guidelines when purchasing hydrogen:

- **Hydrogen Generator**—Because it minimizes the amount of hydrogen present and reduces the degree of hazard, a hydrogen generator (also called an electrolyzer) is the safest way to purchase hydrogen in the quantity used in gas chromatography/mass spectroscopy systems.

However, to minimize the degree of hazard, operate the hydrogen generator only in a non-explosive environment because hydrogen buildup can be ignitable. Thus, your ventilation system for the room or lab hood where the hydrogen generator operates must maintain an air exchange rate at least two orders of magnitude greater than the maximum hydrogen production rate of the hydrogen generator. Follow the manufacturers' directions about proper use and maintenance of the regulator.

To prevent the possibility of releasing hydrogen, set the hydrogen generator to shut down if:

- There is a loss of flow to the ventilation system
- A hydrogen detector alarms at 25% of the lower flammable limit of hydrogen in air.

Vent the oxygen exhausted by the electrolyzer to the outside as well.

- **Hydrogen Cylinder**—Hydrogen can be delivered in standard laboratory gas bottles or cylinders. These cylinders have a limited amount of hydrogen in them and are a safe way to transport and store hydrogen. Always secure, compressed hydrogen gas cylinders, like all compressed gas cylinders, in an upright position, ideally with a non-combustible chain or cable. If the cylinder falls over, the valve can fall off, causing the pressurized cylinder to take off like a rocket, leading to the release of hydrogen and possibly an explosion, severe injury, or death. Never crack a hydrogen cylinder valve to remove dust or dirt from fittings prior to attaching a regulator, as there is a risk of self-ignition.

## Properly Storing Hydrogen

Storing and handling compressed hydrogen gas and cryogenic liquid hydrogen present potential health and safety hazards. Using proper storage and handling techniques is essential to maintaining a safe work environment.

Use the following guidelines when storing hydrogen:

- Store spare hydrogen gas cylinders outside and away from doors, windows, building air intake vents, structures, and vehicle routes. This precaution applies when the hydrogen is or is not in use. Indoor storage of spare hydrogen cylinders has special requirements, which are beyond the scope of this document. Documentation for each vessel should include a description of the vessel, a list of available drawings or other documents, the most recent inspection results, and the responsible person's name.
- Prevent spare cylinders from toppling by wrapping them with chains. The chains should also be protected against corrosion and excessive heat.
- Separate spare hydrogen cylinders from oxidizing gases (such as oxygen) with a 5 ft. (1.5 m) tall fire barrier with a half-hour fire rating or place the cylinders at least 20 ft. (6 m) apart.
- When moving hydrogen cylinders:
  - Remove the regulator and replace the cylinder valve cap before moving.
  - Move cylinders on cylinder carts or with other appropriate transport devices.
  - Never roll or drop a cylinder and never lift a cylinder by its protective cap.
- Bulk hydrogen systems include either gaseous or liquid hydrogen in fixed installations; in some gas systems a semi-permanent trailer (tube trailer) can be used. Storage vessels for compressed hydrogen gas or liquid hydrogen should be designed, constructed, tested, and maintained in accordance with applicable codes and standards. Bulk hydrogen systems represent a level of complexity again which is beyond the scope of this document; however some general guidelines are provided.
- The bulk hydrogen storage system should not be located beneath electric power lines, close to other flammable gases/liquids, or close to public areas. It should be readily accessible to authorized personnel and delivery equipment, but protected from physical damage or tampering.
- As liquid hydrogen systems also have a cryogenic hazard, additional safety considerations for the use of cryogenic liquids might be necessary.

## Hydrogen Safety Codes, Standards and References

The following list of safety codes, standards, and references is in no way an exhaustive list. In fact, there may be federal, state, or local codes that apply to your specific location. Check with all appropriate agencies with jurisdiction before installing or using a hydrogen system.

- Air Products Safetygram #4 Gaseous Hydrogen
- ANSI/AIAA standard for hydrogen safety guidelines is AIAA G-095-2004, Guide to Safety of Hydrogen and Hydrogen Systems
- ASME B31.1, Power Piping Code
- ASME B31.3, Process Piping Code

- ASME B31.8, Gas Transmission and Distribution Systems
- BCGA Code Of Practice CP4 Industrial Gas Cylinder Manifolds and Gas Distribution Pipework
- BCGA Code Of Practice CP33 The Bulk Storage of Gaseous Hydrogen at Users' Premises
- CGA G-5, Hydrogen
- CGA G-5.4, Standard for Hydrogen Piping Systems at Consumer Locations
- CGA G-5.5, Hydrogen Vent Systems
- CGA G-5.6, Hydrogen Pipeline Systems
- CGA G-5.8, High Pressure Hydrogen Piping Systems at Consumer Locations.
- FM Global Property Loss Prevention Data Sheets 7-50: Compressed Gases in Cylinders
- FM Global Property Loss Prevention Data Sheets 7-91: Hydrogen
- IGC Doc 121/04/E, Hydrogen Transportation Pipelines System Design Features
- NASA
- NSS 1740.16 Safety Standard For Hydrogen And Hydrogen Systems Guidelines for Hydrogen System Design, Materials Selection, Operations, Storage, and Transportation
- NFPA 52, Vehicular Fuel Systems Code
- NFPA 55, Standard for the Storage, Use, and Handling of Compressed Gases and Cryogenic Fluids in Portable and Stationary Containers, Cylinders, and Tanks, 2005 Edition
- NFPA 68, Standard on Explosion Protection by Deflagration Venting
- NFPA 70, National Electrical Code
- NFPA 497, Recommended Practice for the Classification of Flammable Liquids, Gases, or Vapors and of Hazardous (Classified) Locations for Electrical Installations in Chemical Process Areas
- NFPA 13, Standard for the Installation of Sprinkler Systems
- NFPA 45, Standard on Fire Protection for Laboratories Using Chemicals
- NFPA 55, Standard for the Storage, Use, and Handling of Compressed Gases and Cryogenic Fluids in Portable and Stationary Containers, Cylinders, and Tanks
- NFPA 68, 2007 Standard on Explosion Protection by Deflagration Venting
- NFPA 69, Standard on Explosion Prevention Systems
- NFPA 91, Standard for Exhaust Systems for Air Conveying of Vapors

- NFPA 255, Standard Method of Test of Surface Burning Characteristics of Building Materials

OSHA 29CFR1910.103 1910.103 Hydrogen

## Hazardous Substances Precautions



**WARNING** Before using hazardous substances (toxic, harmful, and so on), read the hazard indications and information reported in the applicable Material Safety Data Sheet (MSDS.) Use Personal protection according to the safety requirements.

## Venting Toxic Gases

When analyzing toxic compounds be aware that during the normal operation of the GC some of the sample might be vented outside the instrument through the inlet and detector exits; therefore, make sure to vent the exhaust gases to a fume hood. Consult local Environmental and Safety Regulations for instructions in exhausting fumes from your system.

## Liquid Nitrogen Safety Precautions

Liquid nitrogen is a colorless, odorless, extremely cold liquid and gas under pressure. It can cause rapid suffocation when concentrations are sufficient to reduce oxygen levels below 19.5%. A Self Contained Breathing Apparatus (SCBA) might be required. Contact with liquid or cold vapors can cause severe frostbite. Cold vapors in the air will appear as a white fog due to condensation of moisture. Oxygen concentrations must be monitored in the release area. All cryogenic liquids produce large volumes of gas when they vaporize.



**WARNING** Before using Liquid Nitrogen, read the hazard indications and the instructions reported in the Safety sheet supplied by the manufacturer, with reference to the CAS number (Chemical Abstract Service) 7727-37-9.

Use personal protection:

- **Protective gloves:** Loose fitting thermal-insulated or leather gloves.
- **Eye protection:** Full face shield and safety glasses are recommended.
- **Other protective equipment:** Safety shoes when handling containers. Long sleeved shirts and trousers without cuffs. Work clothing that sufficiently prevents skin contact should be worn.

## Carbon Dioxide Safety Precautions

Carbon dioxide is a colorless, cryogenic liquid. At low concentrations, is odorless. At higher concentrations carbon dioxide will have a sharp, acidic odor. At concentrations between 2 and 10%, Carbon dioxide can cause nausea, dizziness, headache, mental confusion, increased blood pressure, and increased respiratory rate. If the gas concentration reaches 10% or more, suffocation and death can occur within minutes. Contact with the cold gas can cause freezing of exposed tissue. Moisture in the air could lead to the formation of carbonic acid that can be irritating to the eyes. All forms of carbon dioxide are noncombustible. Carbon dioxide is heavier than air and should not be allowed to accumulate in low lying areas.

**WARNING** Before using carbon dioxide, read the indications of hazard and the instructions reported in the Safety sheet supplied by the manufacturer with reference to the CAS number (Chemical Abstract Service) 124-38-9.



Use personal protection:

- **Protective gloves:** Loose fitting thermal insulated or leather gloves.
- **Eye protection:** Full face shield and safety glasses are recommended.
- **Other protective equipment:** Safety shoes when handling containers. Long sleeved shirts and trousers without cuffs. Work clothing that sufficiently prevents skin contact should be worn.

# Installation

This chapter provides instructions for installing and connecting the TRACE 1300/TRACE 1310 system.

## Contents

- Positioning the TRACE 1300/TRACE 1310
- Installing the External Accessories
- Installing the Injector and Detector Modules
- Making the Gas Supply Plumbing Connections
- Connecting the Oven Cryogenic System
- Connecting the PTV/PTVBKF Cryogenic System
- Coupling to a Mass Spectrometer
- Installing the Autosampler
- Installing the Data System Software
- Making Power Connections
- Setting the LAN Communication
- Column Installation Requirements
- Installing the Column the First Time



**IMPORTANT** Before beginning the installation process, your laboratory must be in compliance with the guidelines and the requirements described in the *TRACE 1300/TRACE 1310 Preinstallation Requirements Guide*.

When you receive your TRACE 1300/TRACE 1310, an authorized Thermo Fisher Scientific Field Service Engineer (FSE) will configure the system for you. However, if you need to reconfigure your system, repeat the following steps:

- 1. Place the TRACE 1300/TRACE 1310 on the workbench.**
  - a. See [“Positioning the TRACE 1300/TRACE 1310”](#) on [page 3](#).
- 2. Place and connect the external accessories.**
  - a. See [“Installing the External Accessories”](#) on [page 3](#).
- 3. Place and connect the Injector and Detector modules.**
  - a. See [“Installing the Injector and Detector Modules”](#) on [page 4](#).
- 4. Plumb the gas supply and perform test for leaks.**
  - a. See [“Making the Gas Supply Plumbing Connections”](#) on [page 7](#).
  - b. For system test information, refer to the *TRACE 1300 and TRACE 1310 User Guide* and to the *User Guide* of any accessories you have connected to your GC system.
- 5. Connect coolant to cryogenic system (if present).**
  - a. See [“Connecting the Oven Cryogenic System”](#) on [page 13](#).
- 6. Couple a TRACE 1300/TRACE 1310 version for MS to a mass spectrometer.**
  - a. Perform this step to couple a GC MS version with a Thermo Scientific mass spectrometer. See [“Coupling to a Mass Spectrometer”](#) on [page 23](#).
- 7. Install the autosampler (optional) on the GC.**
  - a. See [“Installing the Autosampler”](#) on [page 32](#).
- 8. Install the Data System software.**
  - a. See [“Installing the Data System Software”](#) on [page 42](#).
- 9. Connecting Power to the GC and the External Modules.**
  - a. See [“Making Power Connections”](#) on [page 43](#).
- 10. Set the LAN communication between TRACE 1300/TRACE 1310 and computer.**
  - a. See [“Setting the LAN Communication”](#) on [page 45](#).
- 11. Install the analytical column.**
  - a. Install the column into GC as described in [“Column Installation Requirements”](#) on [page 50](#) and [“Installing the Column the First Time”](#) on [page 57](#).
- 12. Perform column conditioning, leak test, and column evaluation.**
  - a. Perform column conditioning, leak check, and carrier gas flow calibration (column evaluation) following the instruction reported in [“Installing the Column the First Time”](#) on [page 57](#).



## Positioning the TRACE 1300/TRACE 1310

Your laboratory must already be prepared according to the space requirements, and the gas and power supplies must be accessible. Optional equipment should be placed near the TRACE 1300/TRACE 1310 for easier connection.

### ❖ To position the instrument



**WARNING** The main unit of the GC, without injector/detector modules, weighs approximately **35 kg (77 lbs)** when unpacked. Pay attention when lifting the instrument onto the work table.

1. Set the TRACE 1300/TRACE 1310 on the work table. At least **TWO** people should perform this operation, each standing on one side of the instrument and putting their hands near its supporting feet.



**WARNING** The oven vents at the back of the GC discharge hot air up to 450 °C (842 °F) during cooling. Oven exhaust can cause severe burns. Avoid working behind the instrument oven vents during cool-down cycles. There should be at least 30 cm (12 in.) free space behind the instrument to allow the exhaust to dissipate. Do not expose gas tanks or bottles, chemicals, regulators, electrical cables, or other temperature-sensitive items to oven exhaust.

**Note** An optional **oven exhaust kit** can be installed to carry the hot air from the oven vents out of the working area. See [“Installing the Oven Exhaust Kit”](#) on [page 326](#).

2. Look for the yellow label located on the back of the GC under the AC Input Connector. It indicates the power supply (120 Vac or 230 Vac) required by the GC, and must be in compliance with your power source.

**CAUTION** To avoid contact with the hot air from the vents, gather the electrical cables into the cables holder.

## Installing the External Accessories


For additional information, refer to the *TRACE 1300 and TRACE 1310 User Guide*, and to the *User Guides* that come with the accessories you have connected to your GC.

## Installing the Injector and Detector Modules

This section provides the instructions for installing your front/back injector and detector modules into the relevant housing.

The GC is shipped with dummy modules installed into the injector/detector housings. See [Figure 1](#).

**Figure 1.** Injector and Detector Housings and Dummy Modules

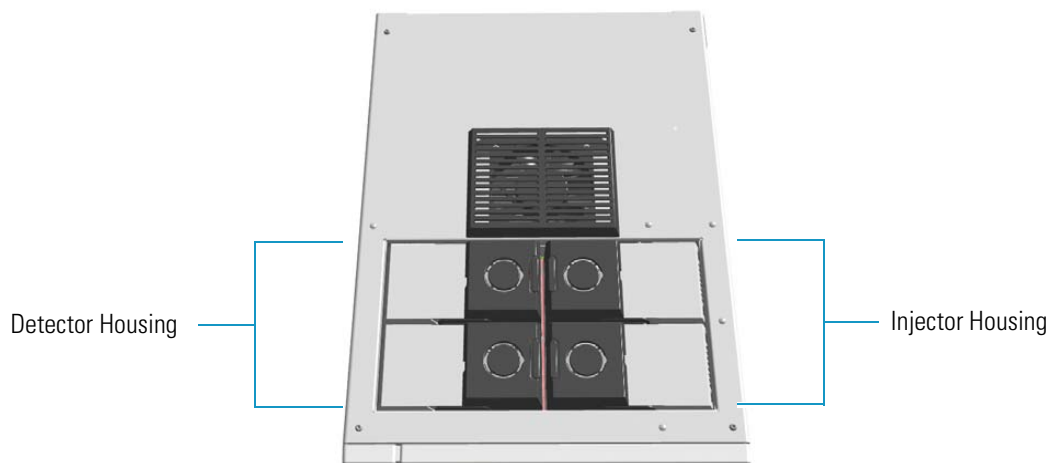

**Note** Where a dummy module is installed, the gas connection is blocked by a plug.

The dummy modules must be removed and replaced by the injector and detector modules as required by the configuration of your TRACE 1300/TRACE 1310.



**CAUTION** Each housing must always be occupied by both modules. If the GC is configured with a single injector/detector module, a dummy module must remain inserted instead of the missing injector/detector module.

The injector and detector modules are shipped with a plug. Before installing the column, remove the plug.

To install the injector and detector modules, see the following procedures:

- [Installing an Injector Module](#)
- [Installing a Detector Module](#)


## Installing an Injector Module

### ❖ To install an injector module

1. Remove the dummy module from the position where the injector module will be installed.

- a. Open the module flap cover.
- b. Using a T20 Torxhead screwdriver, unscrew the two captive fixing screws.
- c. Keeping the dummy module flap cover open, lift up the module from its seat in the injector housing. Place the dummy module on a clean surface.
- d. Remove the gas block plug from the gas connection by unscrewing its fixing screw using a T20 Torxhead screwdriver.

**Figure 2.** Injector Gas Block Plug

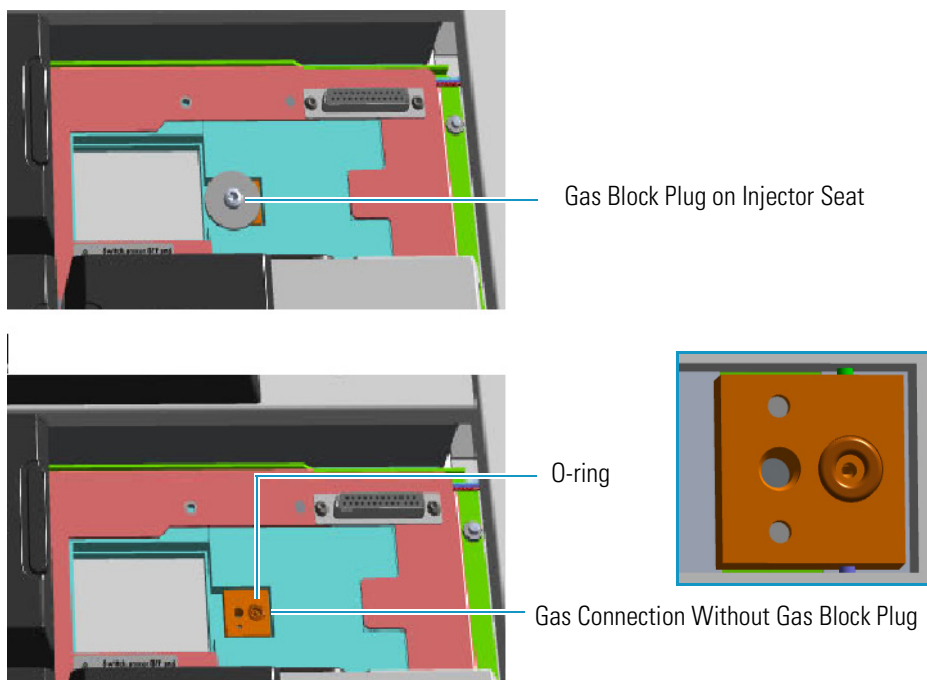

**WARNING** Make sure the o-ring is placed into its seat on the gas connection. See [Figure 2](#). **Do not install the module if the o-ring is missing.**

2. Plug the injector module into the main frame.
  - a. Open the module flap cover.
  - b. Keeping the module flap cover open, place it in its seat. Be sure to insert the 25-pin male connector on the bottom of the module into the 25-pin female connector on the injector seat of the injector housing.
  - c. Use a T20 Torxhead screwdriver to tighten the three captive fixing screws evenly and carefully without overtightening.



**CAUTION** To maintain the correct alignment the screws must be tightened in turn, and each screw must be tightened only a small amount before moving to the next screw. Repeat until all are secure.

## 1 Installation

### Installing the Injector and Detector Modules

- d. Close the module flap cover.
3. Continue the installation following the instructions reporting in the section [“Making the Gas Supply Plumbing Connections”](#) on [page 7](#).

## Installing a Detector Module




**CAUTION** If you are installing a **NPD detector module**, the installation of the NPD Thermionic Source Power Module is required. For the installation details see the section [“Adding a NPD Detector Module”](#) on [page 435](#).

**CAUTION** If you are installing a **Generic Detector Interface** see the section [“Adding a Generic Detector Interface”](#) on [page 472](#).

### ❖ To install a detector module

1. Remove the dummy module from the position where the detector module will be installed.
  - a. Open the module flap cover.
  - b. Using a T20 Torxhead screwdriver, unscrew the two captive fixing screws.
  - c. Keeping the dummy module flap cover open, lift up the module from its seat in the injector/detector housing. Place the dummy module on a clean surface.
  - d. Remove the gas block plug from the gas connections by unscrewing its fixing screw using a T20 Torxhead screwdriver.

**Figure 3.** Detector Gas Block Plug

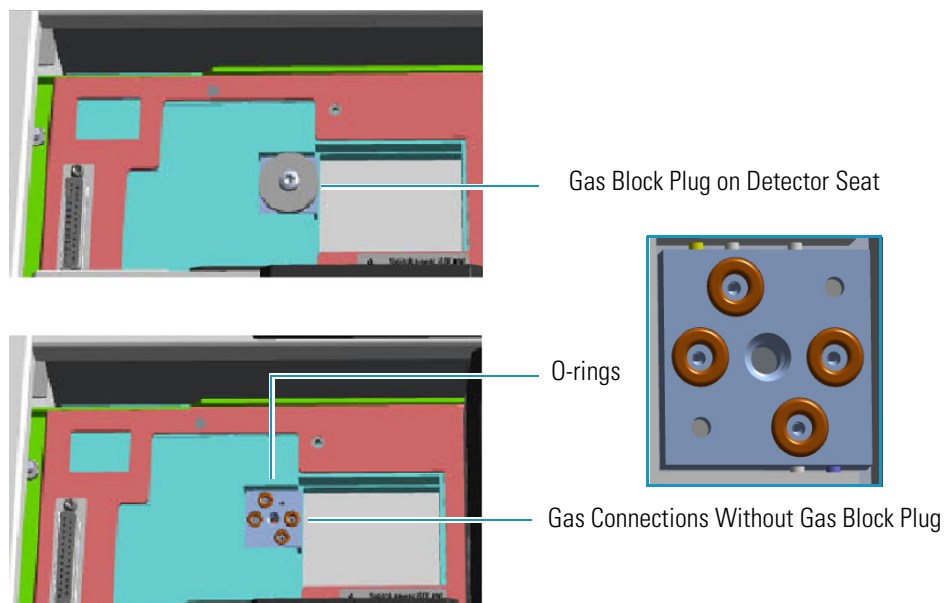



**WARNING** Make sure all the four o-rings are placed into their seats on the gas connection. See [Figure 3](#). **Do not install the module if the o-rings are missing.**

2. Plug the detector module into the main frame.
  - a. Open the module flap cover.
  - b. Keeping the module flap cover open, place the module in its seat.  
Be sure to insert the 25-pin male connector on the bottom of the module into the 25-pin female connector on the detector seat of the detector housing.
  - c. Use a T20 Torxhead screwdriver to tighten the three captive fixing screws evenly and carefully without overtightening.



**CAUTION** To maintain the correct alignment the screws must be tightened in turn, and each screw must be tightened only a small amount before moving to the next screw. Repeat until all are secure.

- d. Close the module flap cover.
3. Continue the installation following the instructions in the section [“Making the Gas Supply Plumbing Connections”](#) on [page 7](#).

## Making the Gas Supply Plumbing Connections

Make plumbing connections between the gas supply lines and the GC gas inlets on the instrument back using the proper fittings



**CAUTION** DO NOT loosen or remove caps from the TRACE 1300/TRACE 1310 until you have purged gas lines, and you are ready to connect them. Loosening or removing caps early will contaminate instruments and filters.

**WARNING** Before using gases, carefully read the hazard indications and information reported in the Safety Sheet supplied by the manufacturer referring to the CAS (Chemical Abstract Service) number. It is the user's responsibility to see that all local safety regulations for the use of gases are obeyed.



All Thermo Fisher Scientific gas chromatographs use an inert gas as the carrier gas. If you wish to use hydrogen as a carrier gas, you must install a hydrogen sensor. Contact a Thermo Fisher Scientific sales representative if you plan to use hydrogen as the carrier gas in your new TRACE 1300/TRACE 1310. If you don't have the hydrogen sensor, you **must** use an inert carrier gas.

Thermo Fisher Scientific FSEs are not authorized to install or repair any instrument using hydrogen as a carrier gas unless the instrument is equipped with the appropriate sensor.


## 1 Installation

### Making the Gas Supply Plumbing Connections

#### ❖ To plumb the gas supply


1. By now you must have done the following:
  - a. Completely built your gas supply lines including any traps, tees, and extra tubing to allow about 40 cm (16 in.) of slack in the line.
  - b. Purged the gas line after every tube cut to remove any debris or contaminants.
  - c. Ensured the gas supply is turned off.
2. Connect the gas lines.
  - a. TRACE 1300/TRACE 1310 is provided with six gas inlet ports for the connection of carrier and detector gases. See [Figure 4](#).

**Figure 4.** Gases Inlet Connections

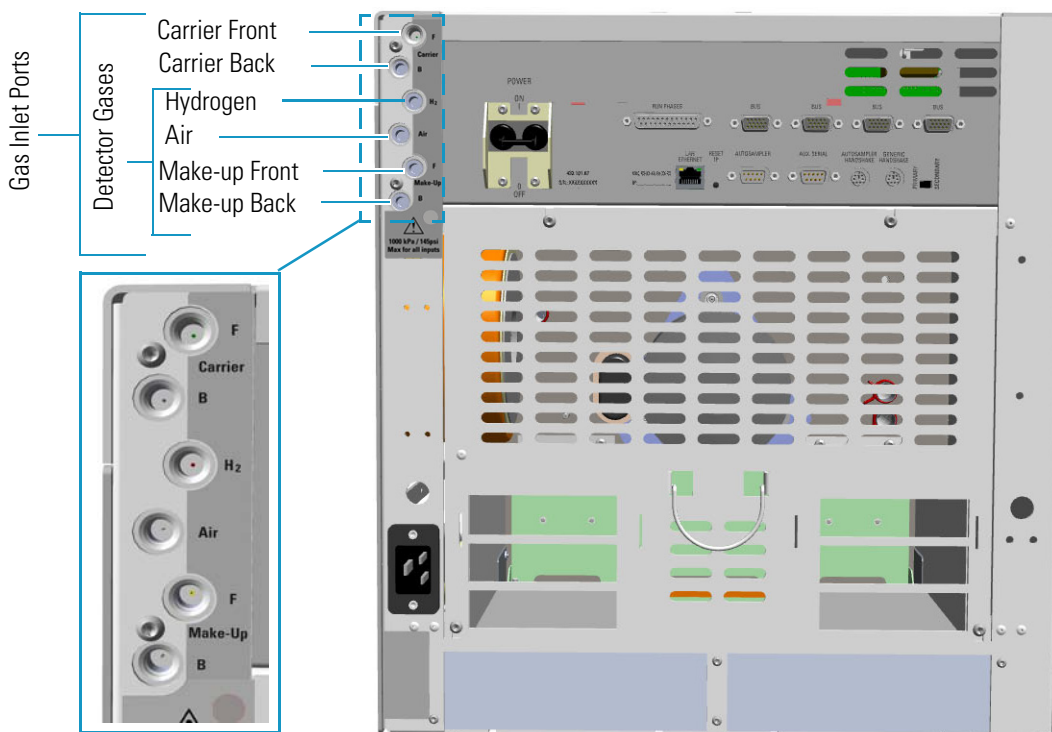

**Note** Use the 1/8-in. Swagelok fittings provided on the gas inlet ports (see [Figure 5](#)) to connect the gas lines.

**Figure 5.** Fittings for Gas Inlets Connection

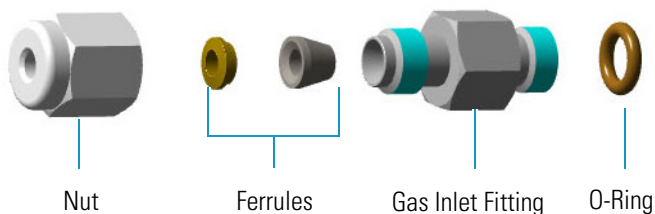




**CAUTION** Inside each 1/8-in. inlet manifold is a fritted filter. To keep the gas line from touching, and possibly damaging the filter, extend the tubing only 5-mm past the front ferrule. This ensures the tubing does not touch the filter.

- b. If present, remove the cap from the gas inlet port to connect. See [Figure 6](#).

**Figure 6.** Cap Removal

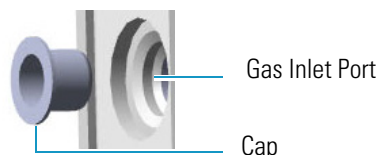


- c. If not already installed, screw the gas inlet fitting provided into the gas inlet port interposing the o-ring. See [Figure 7](#).

**Figure 7.** Gas Inlet Fitting Installation




- d. Connect the gas line to the front/back carrier gas inlet. Use a 7/16-in. wrench for tightening the fittings. See [Figure 8](#).

**Figure 8.** Carrier Gas Line Connection

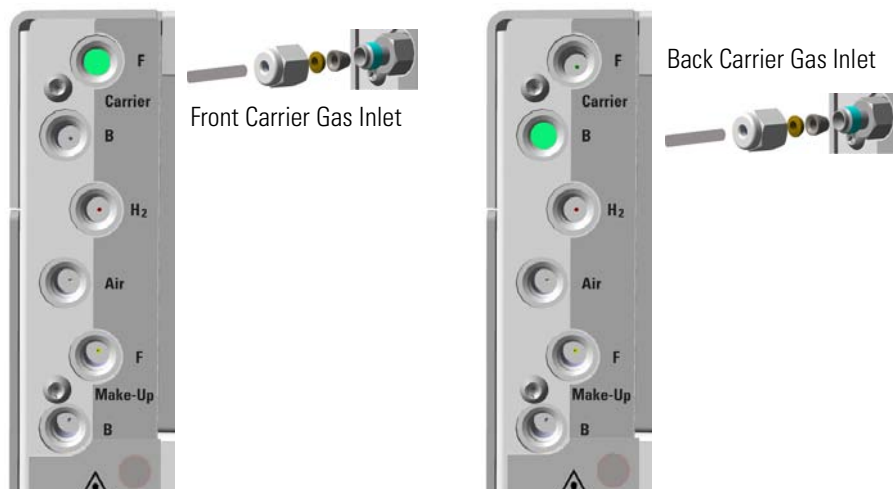


- e. Connect the gas line to the front/back detector gas inlet. Use a 7/16-in. wrench for tightening the fittings. See [Figure 9](#), [Figure 10](#), [Figure 11](#), and [Figure 12](#) to properly connect each detector.


**Figure 9.** FID/NPD Detector Gas Lines Connection

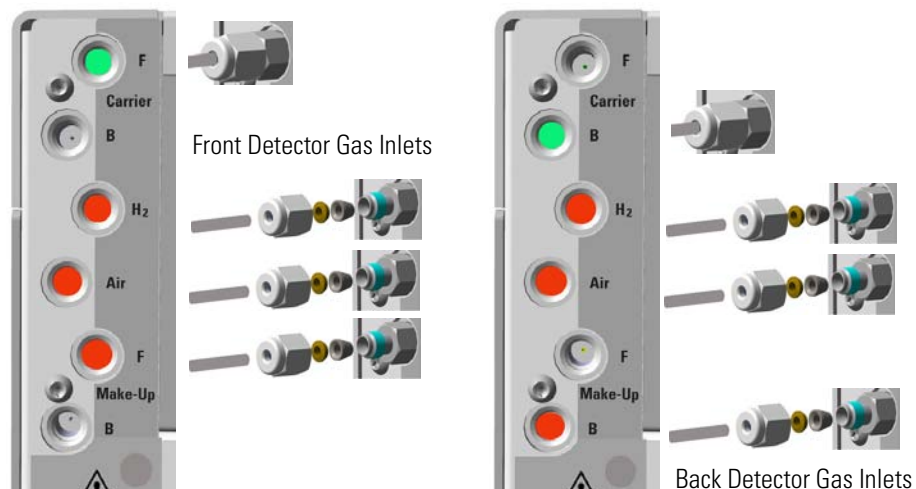

**Figure 10.** ECD Detector Gas Lines Connection

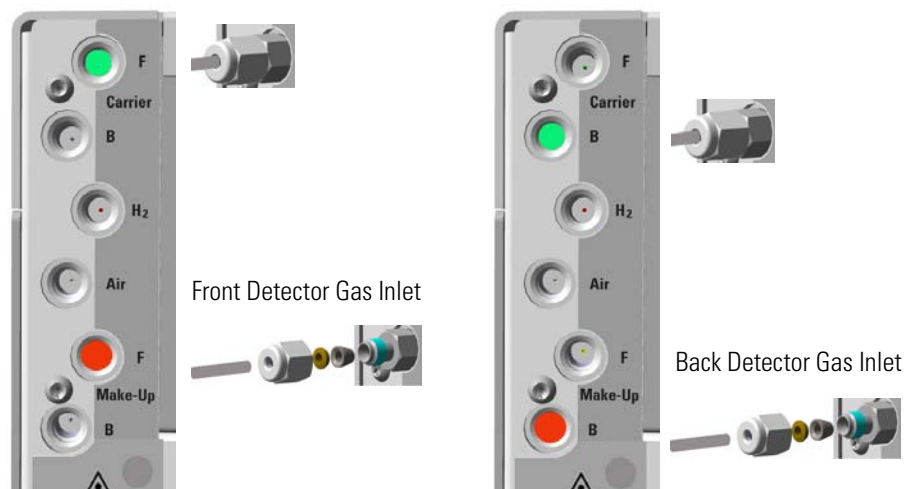


**Figure 11.** TCD Detector Gas Lines Connection

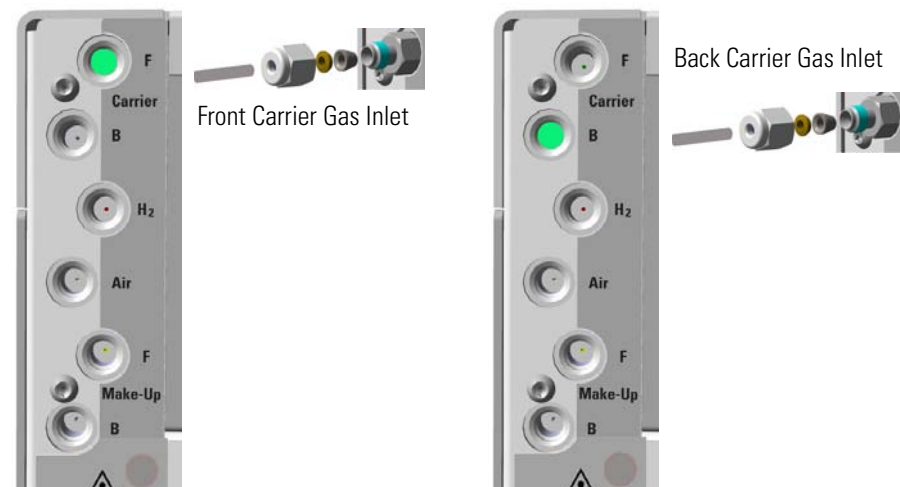



**Figure 12.** FPD Detector Gas Lines Connection

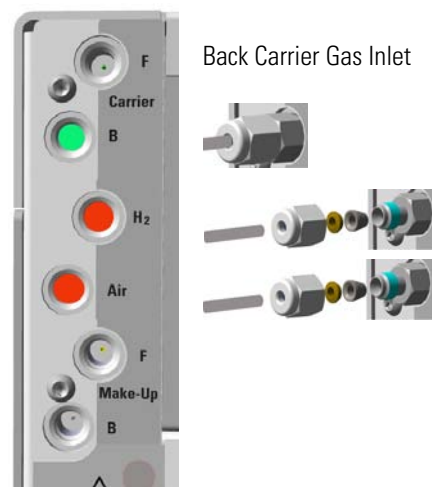

- f. Be sure to complete all inlet and detector connections before turning on the gas supply.

3. Supply the gas lines.



**IMPORTANT** The maximum nominal inlet pressure for all the inputs is 1000 kPa (145 psig), as indicated on the label under the gas inlets ports on the back of the GC. The working inlet pressure range is from 400 kPa (58 psig) to 1000 kPa (145 psig).

**Note** If auxiliary carrier gases are required, the installation of the **Auxiliary Gas** module is required. For the installation details, see the section [“Adding an Auxiliary Gas System”](#) on [page 521](#).

## Testing for Leaks

Once you have connected the gas supplies to the GC, you need to test the gas supply lines for leaks.



**CAUTION** Before performing a leak test make sure that the GC is powered off and the power cable is unplugged from the AC Input connector (Mains socket) and from the wall outlet.

### ❖ To perform a leak test

1. Be sure that GC is powered off.
2. Open and set the gas supply.
  - a. Open the gas supply.
  - b. Set the carrier gas pressure to approximately 50 kPa (7 psi) higher than the maximum pressure of the GC regulator.
  - c. Set the detector gas pressures to approximately 1000 kPa (145 psig) if available, otherwise set the maximum pressure allowed, for example 500 kPa (72.5 psi).
3. Check for leaks.
  - a. Use a handheld electronic leak detector (Thermo Scientific GLD Pro leak detector or equivalent) to check each fitting for leaks.
  - b. If you detect a leak, tighten the connection and retest it.
  - c. Repeat this process until all connections are leak free.


## Connecting the Oven Cryogenic System

If the oven of your TRACE 1300/TRACE 1310 is equipped with a cryogenic cooling option, you need a supply of coolant, such as carbon dioxide (CO<sub>2</sub>) or liquid nitrogen (N<sub>2</sub>).

### Oven Cryo System with Carbon Dioxide Connection

This Oven cryo system consists of a tube for carbon dioxide inserted into the oven and fitted with a 1/8-in. solenoid valve for carbon dioxide mounted into the back of the GC. Carbon dioxide must be supplied in a high-pressure cylinder with a dip tube. It is your responsibility to ensure the delivery connection from the carbon dioxide cryogenic supply is adaptable to 1/8-in. tubing. [Figure 13](#) shows the proper carbon dioxide container configuration.

**Figure 13.** Carbon Dioxide Container

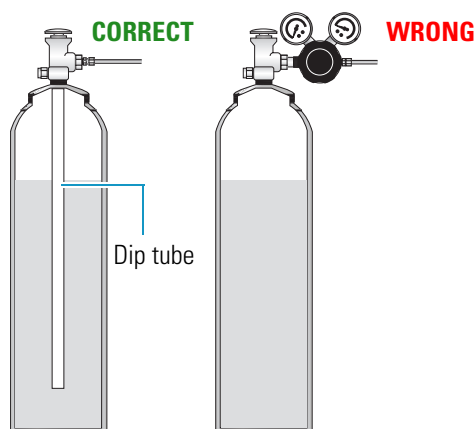


**CAUTION** High pressures and extremely low temperatures make pressurized carbon dioxide a hazardous material. High concentrations of carbon dioxide are dangerous. To avoid injury, always follow the safety precautions and delivery system design recommended by your gas supplier.

#### ❖ To connect the cryogenic system to the carbon dioxide supply

1. Connect the proper end of the cryo supply tube to the 1/8-in. NPT connection of the solenoid valve using the proper nut and ferrule. Use a 7/16-in. wrench to tighten the fittings. See [Figure 14](#).

**Figure 14.** Oven Cryo Supply Tube for Carbon Dioxide Connection

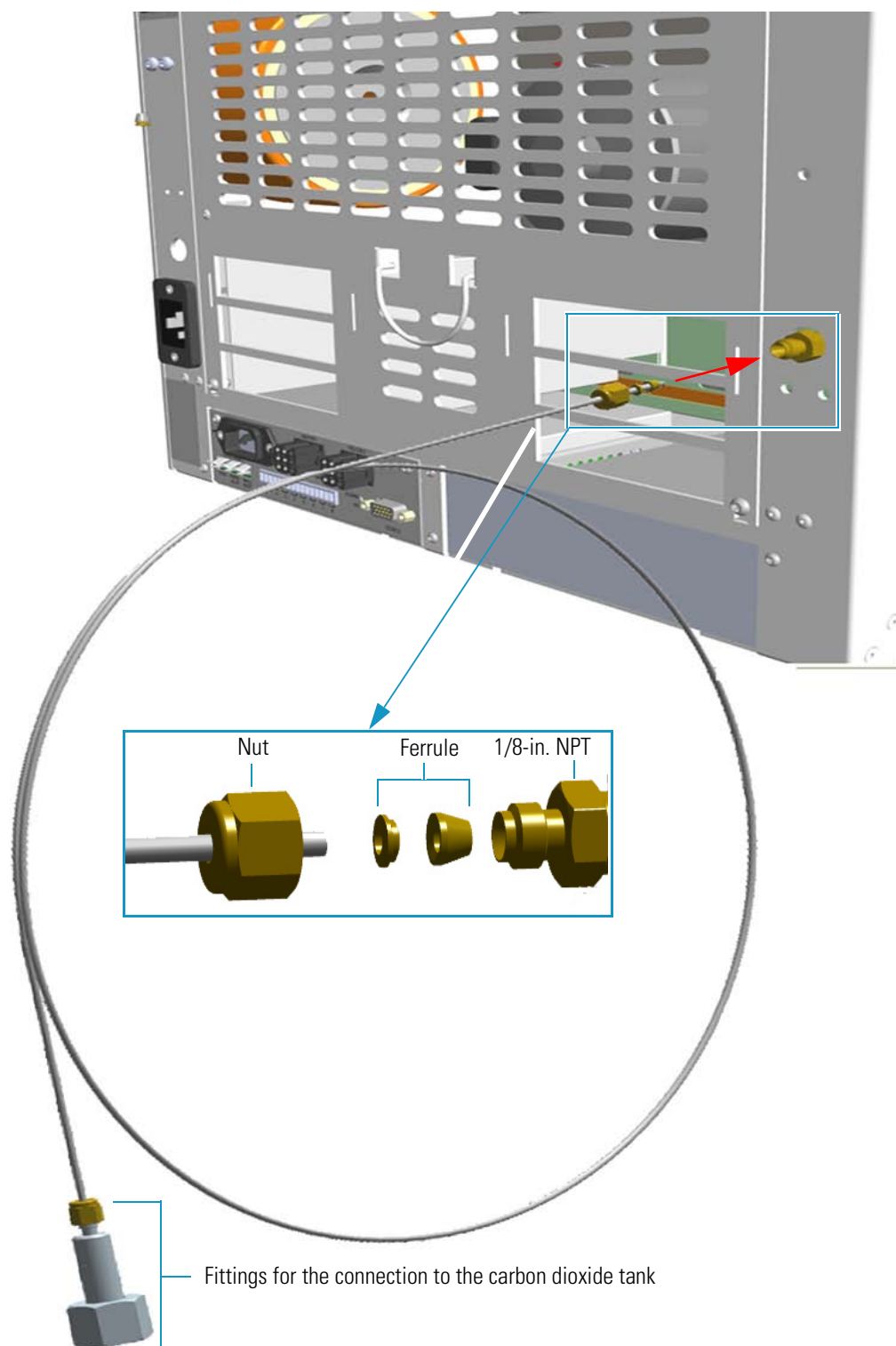

2. Connect the other end of the cryo supply tube to the coolant container using the appropriate nut and ferrule.


3. Connect the cryo valve to the **Aux Temperature/Cryo Module**.



**IMPORTANT** The Auxiliary Temperature/Cryo must be powered at the same line voltage of the main GC system.

- a. Connect the cryo solenoid valve to the 2-pin connector marked **Cryo Valves - Oven** using the cable provided. See [Figure 15](#).

**Figure 15.** Cryo Valves: Oven

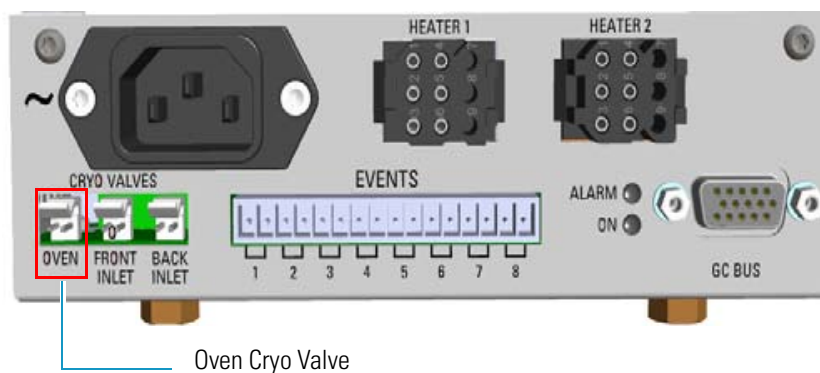

- b. Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.
- c. Plug the power cable to the AC Input connector on the front of the module and to the wall outlet. The LED marked **On** lights up after the GC power on.

**Note** For further details regarding the installation of the Aux Temperature/Cryo module, see the section [“Adding an Aux Temperature/Cryo Module”](#) on [page 443](#).

## Oven Cryo System with Liquid Nitrogen Connection

This Oven cryo system consists of a tube for liquid nitrogen inserted into the oven and fitted with a 1/8-in. solenoid valve for liquid nitrogen mounted into the back of the GC.

The tube for liquid nitrogen is covered with an insulating material tube.

Liquid Nitrogen must be supplied at a pressure of 1.5 bar (150 kPa; 21.75 psig). Plumbing to the GC should be 1/4-in. copper or stainless steel tubing with insulation. It is your responsibility to ensure the delivery connection from the liquid nitrogen cryogenic supply is adaptable to 1/4-in. tubing. The liquid nitrogen cryogenic valve on the GC is a 1/4-in. Swagelok fitting.

Figure 24 shows the proper configuration for a liquid nitrogen tank.

**Figure 16.** Liquid Nitrogen Tank Configuration




**WARNING** High pressures and extremely low temperatures make liquid nitrogen a hazardous material. High concentrations of liquid nitrogen in the air can be an asphyxiation hazard. To avoid injury, always follow the safety precautions and delivery system design recommended by your gas supplier.

### ❖ To connect the cryogenic system to the liquid nitrogen supply

1. Connect the proper end of the cryo supply tube to the 1/8-in. NPT connection of the solenoid valve using the proper nut and ferrule. Use a 7/16-in. wrench to tighten the fittings. See Figure 17.

**Figure 17.** Oven Cryo Supply Tube for Liquid Nitrogen Connection

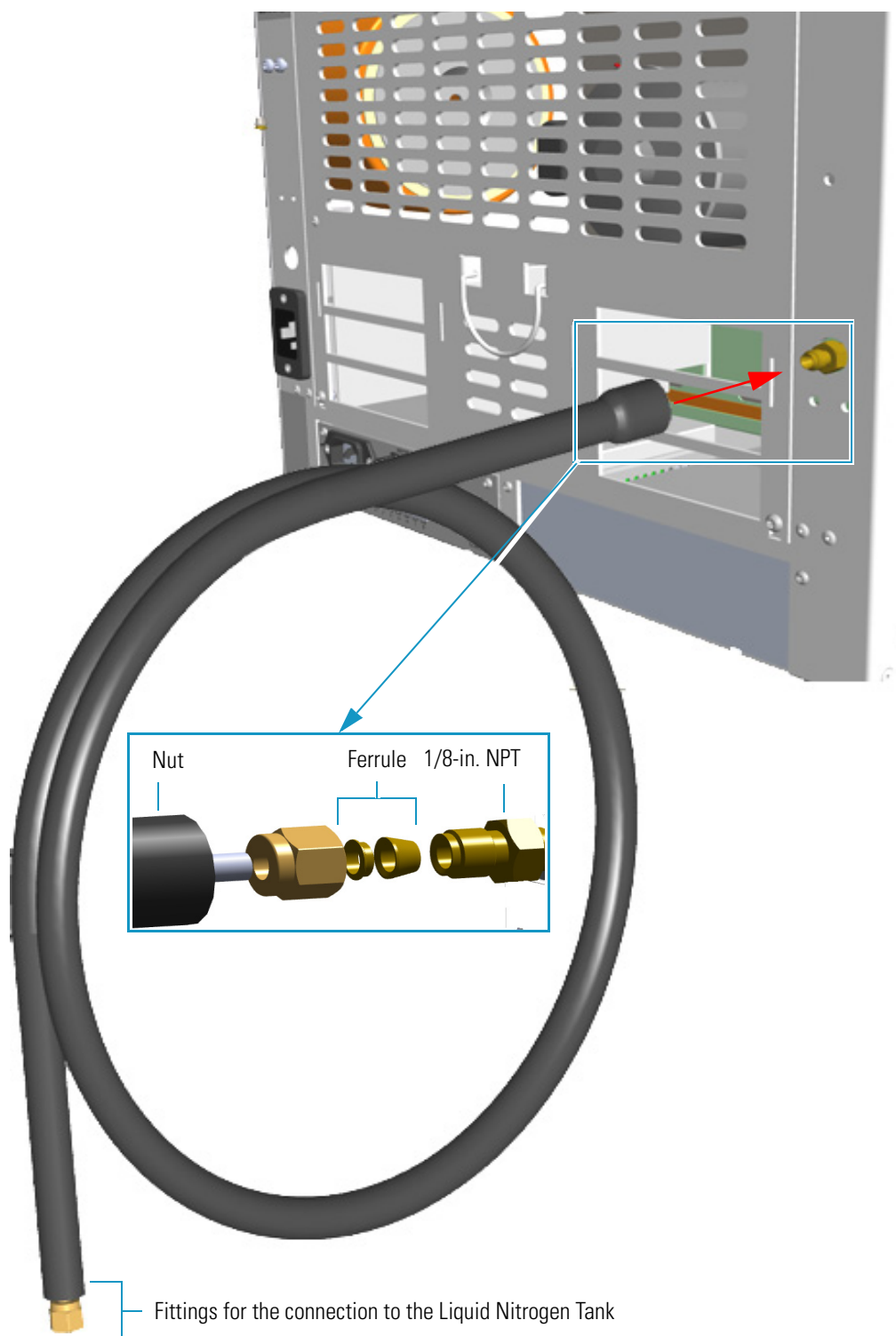

2. Connect the other end of the cryo supply tube to the coolant container using the appropriate nuts and ferrules.

## 1 Installation

### Connecting the PTV/PTVBKF Cryogenic System


3. Connect the cryo valve to the **Aux Temperature/Cryo Module**.



**IMPORTANT** The Auxiliary Temperature/Cryo must be powered at the same line voltage of the main GC system.

- a. Connect the cryo solenoid valve to the 2-pin connector marked **Cryo Valves - Oven** using the cable provided. See [Figure 18](#).

**Figure 18.** Cryo Valves: Oven

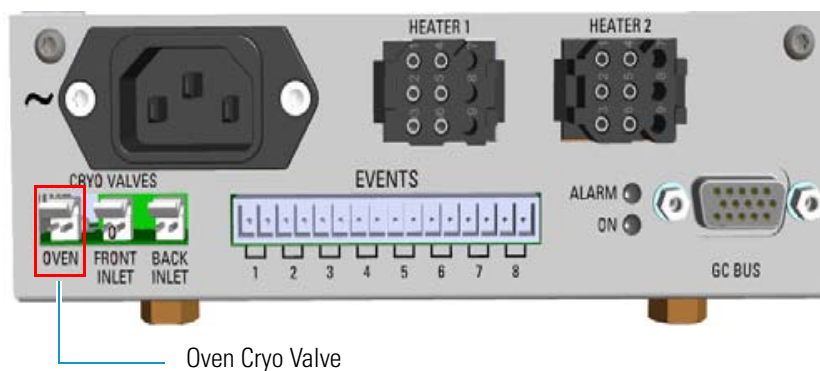

- b. Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.
- c. Plug the power cable to the AC Input connector on the front of the module and to the wall outlet. The LED marked **On** lights up after the GC power on.

**Note** For further details regarding the installation of the Aux Temperature/Cryo module, see the section [“Adding an Aux Temperature/Cryo Module”](#) on [page 443](#).

## Connecting the PTV/PTVBKF Cryogenic System

If the PTV or PTVBKF injector of your TRACE 1300/TRACE 1310 is equipped with a cryogenic cooling option, you need a supply of coolant, carbon dioxide (CO<sub>2</sub>) or liquid nitrogen (LN<sub>2</sub>).

### PTV/PTVBKF Cryo System with Carbon Dioxide Connection

This cryo system consists of a tube for Carbon Dioxide inserted into the PTV/PTVBKF module, and fitted with a 1/8-in. solenoid valve for carbon dioxide mounted on a bracket fixed on the back of the GC. See [Figure 19](#).




**Figure 19.** Front and Back PTV Cryo System with CO<sub>2</sub>



Carbon dioxide must be supplied by a high-pressure cylinder with a dip tube. It is your responsibility to ensure the delivery connection from the carbon dioxide cryogenic supply is adaptable to 1/8-in. tubing. [Figure 20](#) shows the proper carbon dioxide container configuration.

**Figure 20.** Carbon Dioxide Container

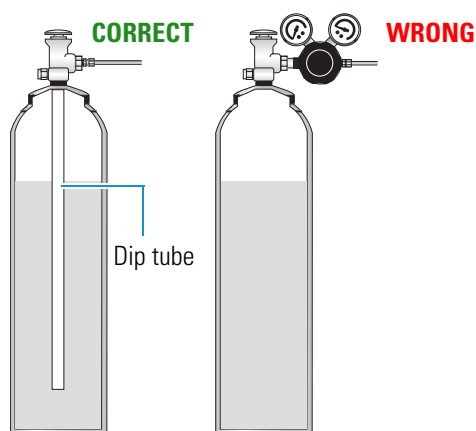

**CAUTION** High pressures and extremely low temperatures make pressurized CO<sub>2</sub> a hazardous material. High concentrations of CO<sub>2</sub> are dangerous. To avoid injury, always follow the safety precautions and delivery system design recommended by your gas supplier.


❖ **To connect the cryogenic system to the carbon dioxide supply**

1. Connect the proper end of the cryo supply tube to the 1/8-in. NPT connection of the solenoid valve using the proper nut and ferrule. Use a 7/16-in. wrench to tighten the fittings. See [Figure 21](#).

## 1 Installation

Connecting the PTV/PTVBKF Cryogenic System

**Figure 21.** PTV/PTVBKF Cryo Supply Tube for Carbon Dioxide Connection

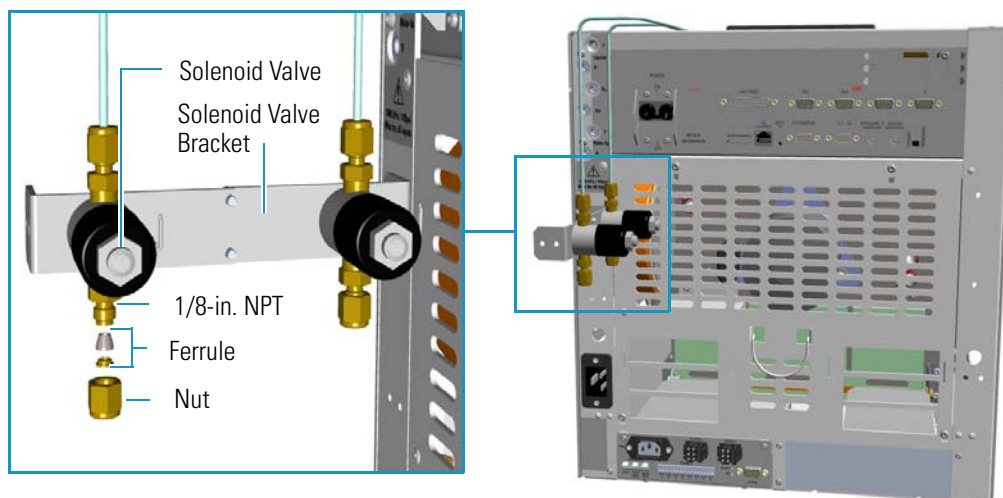


2. Connect the other end of the cryo supply tube to the coolant container using the appropriate nuts and ferrules.
3. Connect the cryo valve to the **Aux Temperature/Cryo Module**.



**IMPORTANT** The Aux Temperature/Cryo module must be powered at the same line voltage of the main GC system.

- a. Connect the 2-pin connector marked **Cryo Valves - Front Inlet** or **Cryo Valves-Back Inlet** to the cryo solenoid valve using the cable provided. See [Figure 22](#).

**Figure 22.** Cryo Valves: Front/Back Inlet

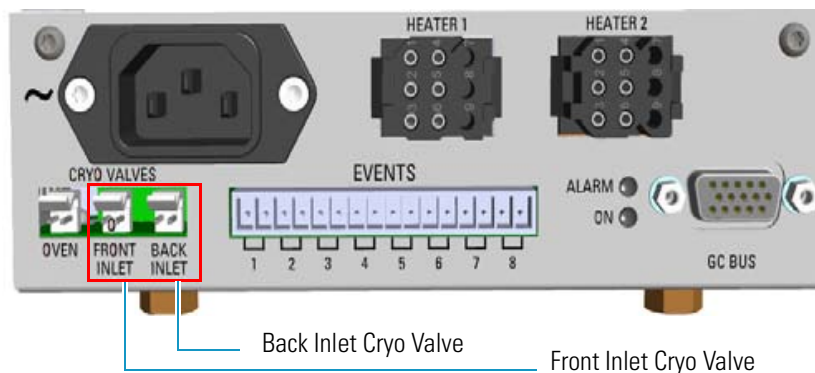


- b. Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.
- c. Plug the power cable to the AC Input connector on the front of the module and to the wall outlet. The LED marked **On** lights up after the GC power on.

**Note** For further details regarding the installation of the Aux Temperature/Cryo module, see the section [“Adding an Aux Temperature/Cryo Module”](#) on [page 443](#).

## PTV/PTVBKF Cryo System with Liquid Nitrogen Connection

This cryo system consists of a tube for liquid nitrogen inserted into the PTV/PTVBKF module and fitted with a 1/8-in. solenoid valve for liquid nitrogen mounted on a bracket fixed on the back of the GC. The tube for liquid nitrogen is covered with an insulating material tube. See [Figure 23](#).

**Figure 23.** Front and Back PTV Cryo System with Liquid Nitrogen

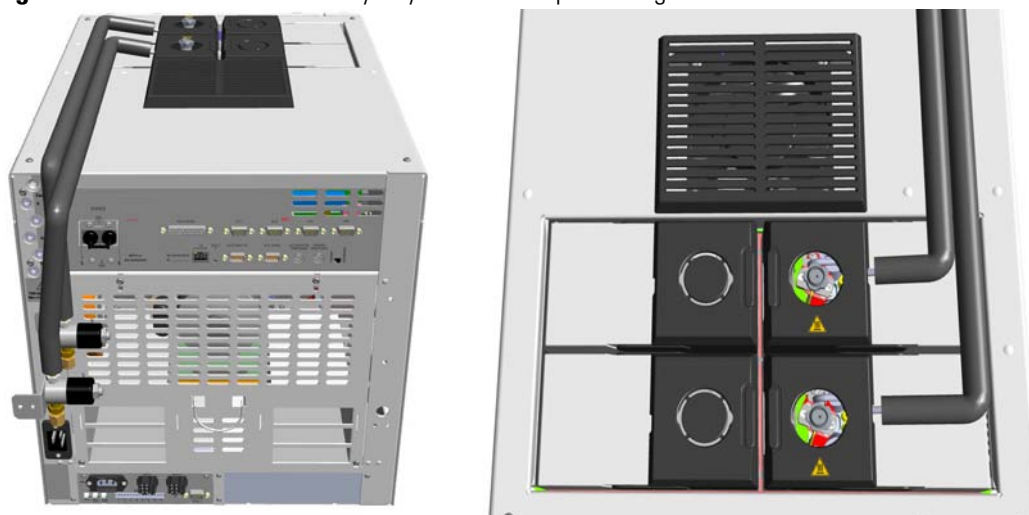

Liquid nitrogen must be supplied at a pressure below 400 kPa (60 psi). Plumbing to the GC should be 1/4-in. copper or stainless steel tubing with insulation.

It is your responsibility to ensure the delivery connection from the Liquid Nitrogen cryogenic supply is adaptable to 1/4-in. tubing. The liquid nitrogen cryogenic valve on the GC is a 1/4-in. Swagelok fitting.

[Figure 24](#) shows the proper configuration for a liquid nitrogen tank.

**Figure 24.** Liquid Nitrogen Tank Configuration



**WARNING** High pressures and extremely low temperatures make Liquid Nitrogen a hazardous material. High concentrations of liquid nitrogen in the air can be an asphyxiation hazard. To avoid injury, always follow the safety precautions and delivery system design recommended by your gas supplier.


## 1 Installation

Connecting the PTV/PTVBKF Cryogenic System

### ❖ To connect the cryogenic system to the liquid nitrogen supply

1. Connect the proper end of the cryo supply tube to the 1/8-in. NPT connection of the solenoid valve using the proper nut and ferrule. Use a 7/16-in. wrench to tighten the fittings. See [Figure 21](#).

**Figure 25.** PTV/PTVBKF Cryo Supply Tube for LN<sub>2</sub> Connection

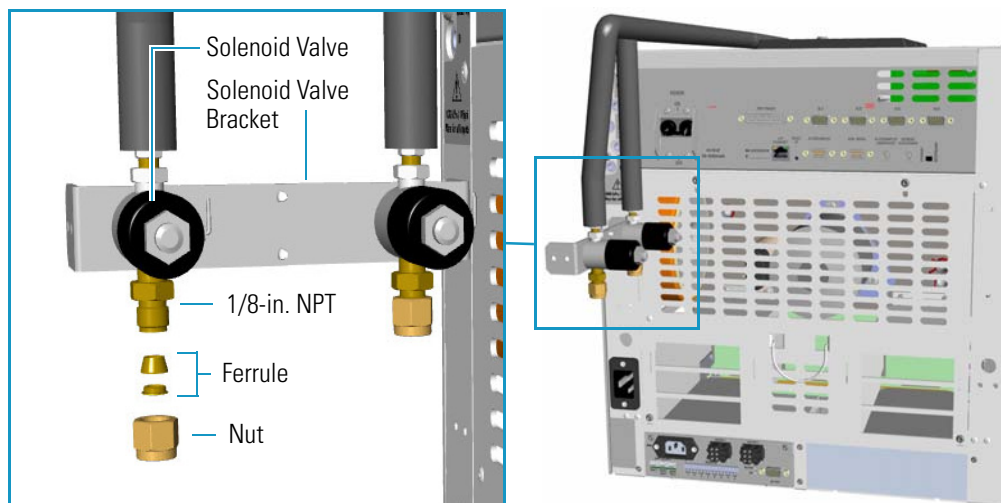


2. Connect the other end of the cryo supply tube to the coolant container using the appropriate nuts and ferrules.
3. Connect the cryo valve to the **Aux Temperature/Cryo Module**.



**IMPORTANT** The Auxiliary Temperature/Cryo must be powered at the same line voltage of the main GC system.

- a. Connect the 2-pin connector marked **Cryo Valves - Front Inlet** or **Cryo Valves-Back Inlet** to the cryo solenoid valve using the cable provided. See [Figure 26](#).

**Figure 26.** Cryo Valves: Front/Back Inlet

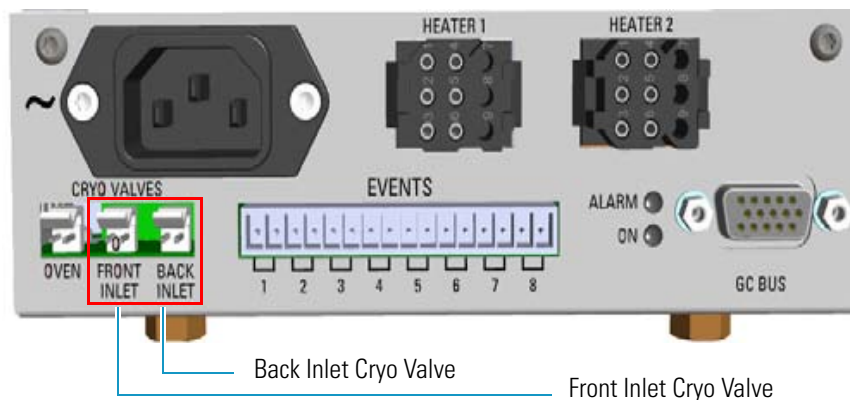

- b. Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.

- c. Plug the power cable to the AC Input connector on the front of the module and to the wall outlet. The LED marked **On** lights up after the GC power on.

**Note** For further details regarding the installation of the Aux Temperature/Cryo module, see the section “Adding an Aux Temperature/Cryo Module” on page 443.

## Coupling to a Mass Spectrometer

This section provides instruction for opening the duct to introduce the transfer line inner tube into the oven of your TRACE 1300/TRACE 1310 version for MS.

According to the mass spectrometer unit to couple, see the following sections:


- “Making the Duct to Couple with an ISQ Series, TSQ 8000 Series, DSQ II, or ITQ Mass Spectrometer” on page 23.
- “Making the Duct to Couple with a TSQ Quantum Mass Spectrometer” on page 26.
- “Making the Duct to Couple with a DFS, IRMS, or ICP-MS Mass Spectrometer” on page 29.

### Making the Duct to Couple with an ISQ Series, TSQ 8000 Series, DSQ II, or ITQ Mass Spectrometer

#### ❖ To make the duct for the transfer line


1. Remove the left side panel.
  - a. Open the front door of the GC.
  - b. Use a T20 Torxhead screwdriver to loosen the left side panel screw from the interior front panel. See Figure 27. Save the screw because it will be reused later.

**Figure 27.** Left Side Panel Fixing Screw

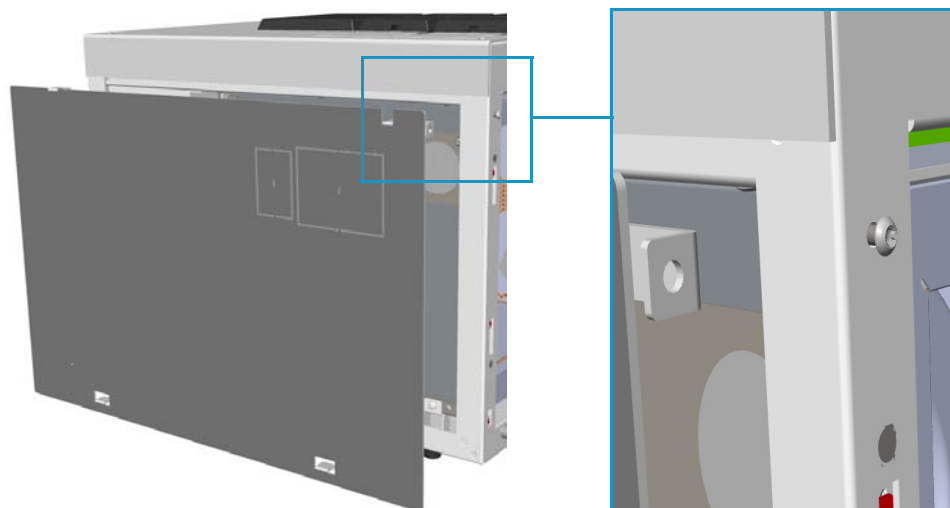


- c. Slide the panel towards the back of the instrument up to the stop.
  - d. Remove the panel by pulling it outward. Be aware that the ground wire is attached to the panel.
2. Prepare the duct for the ISQ Series, TSQ 8000 Series, DSQ II, or ITQ transfer line inner tube.
- a. Remove the partial cut plate on the left panel for that corresponds to your MS. See [Figure 28](#).

**Figure 28.** Left Panel for MS Partial Cut Plate

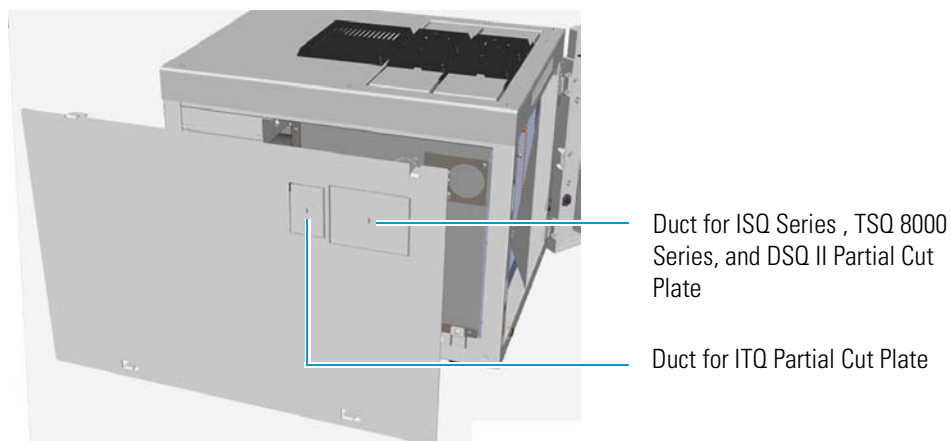

- b. On the exterior wall of the oven box, remove the pre-shaped plug of insulating material from the duct provided. See [Figure 29](#).

**Figure 29.** Pre-shaped Plug Removal

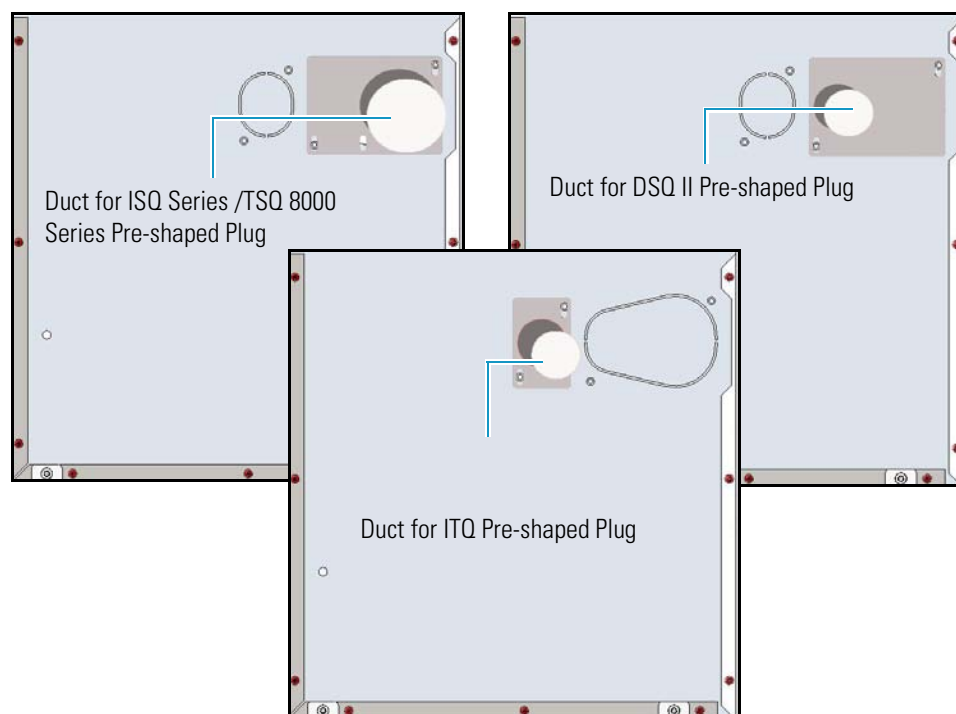




**IMPORTANT** Save the pre-shaped plug of insulating material in a safe place because it could be reused.

- c. On the interior of the oven box remove the partial cut plate from the corresponding duct. See [Figure 30](#).

**Figure 30.** Ducts for the Coupling with ISQ Series, TSQ 8000 Series, DSQ II and ITQ Mass Spectrometer

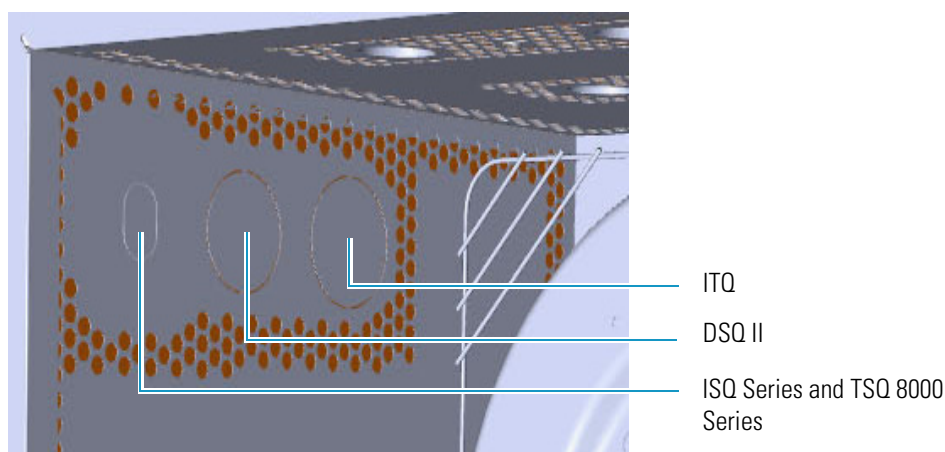

- d. Place the left panel for your MS and attach the screw holding it in place.
3. Introduce the transfer line inner tube into the oven through the duct provided.
  4. Attach the transfer line to the GC column using the proper nut and ferrule.



**CAUTION - INSTRUMENT DAMAGE:** Make sure that the GC column has been conditioned before connecting it to the transfer line. The material released from the column, (column bleed), during conditioning may contaminate the detector.


5. Close the front door of the GC.
6. For DSQ II and ITQ, connect the transfer line heater to the **Aux Temperature/Cryo Module**. See [Figure 31](#).



**CAUTION** An extension cable transfer line is needed when connecting an ITQ mass spectrometer to the Aux Temperature/Cryo module.



**Figure 31.** Transfer Line MS Heaters

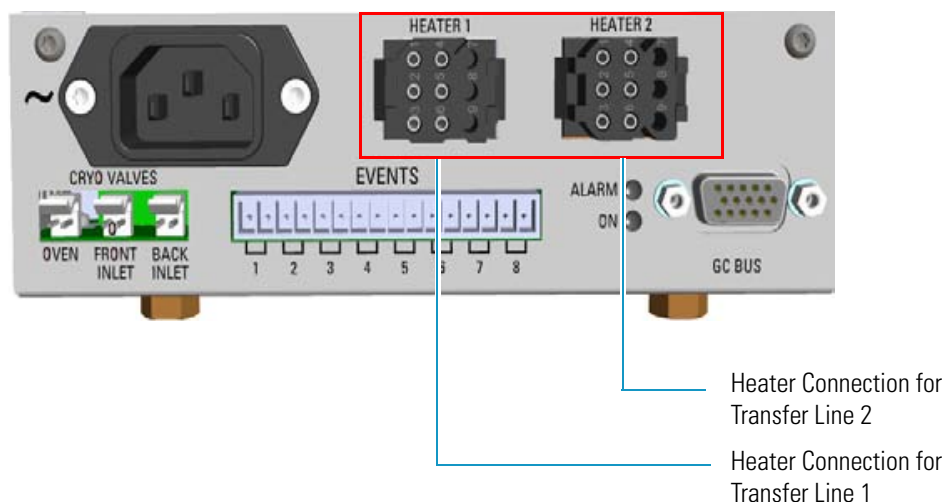


7. To tune and set the ISQ Series, TSQ 8000 Series, DSQ II, and ITQ working conditions, refer to the relevant *User Guide* and *Hardware Manual*.

## Making the Duct to Couple with a TSQ Quantum Mass Spectrometer

### ❖ To make the duct for the transfer line

1. Remove the right side panel.
  - a. Open the front door of the GC.
  - b. Use a T20 Torxhead screwdriver to loosen the right side panel screw from the interior front panel. See [Figure 32](#). Save the screw because it will be reused later.

**Figure 32.** Right Side Panel Fixing Screw

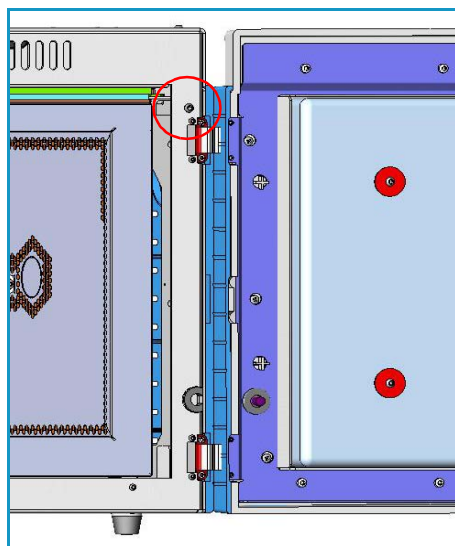


- c. Slide the panel towards the back of the instrument up to the stop.




- d. Remove the panel by pulling it outward. Be aware that the ground wire is attached to the panel.
2. Prepare the duct for the TSQ Quantum transfer line inner tube.
- a. Use a T20 Torxhead screwdriver to loosen the duct plate screws from the GC right panel of the MS. See [Figure 33](#).

**Figure 33.** Right Panel for MS

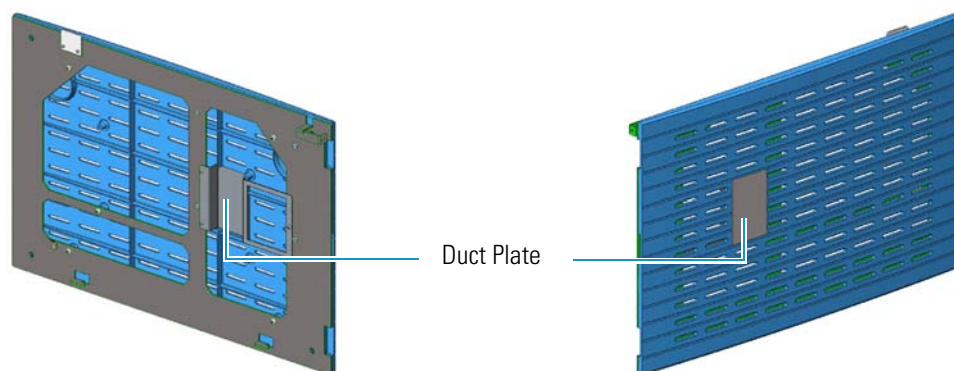

- b. On the exterior wall of the oven box, remove the pre-shaped plug of insulating material from the duct provided. See [Figure 34](#).


**Figure 34.** Pre-shaped Plug Removal

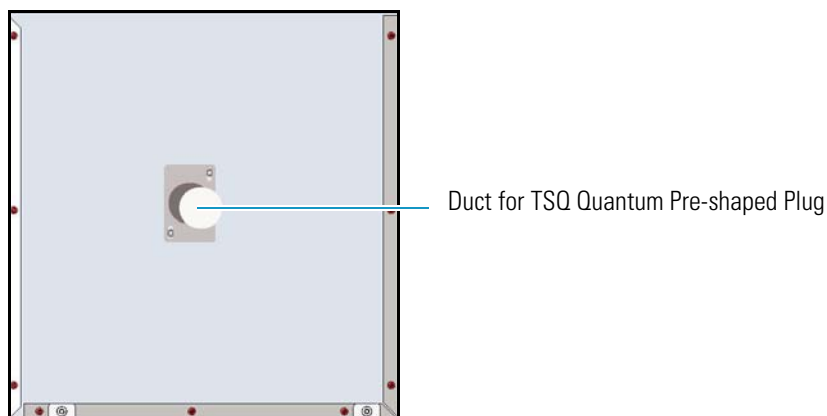

**IMPORTANT** Save the pre-shaped plug of insulating material in a safe place because it could be reused.

- c. On the interior of the oven box remove the partial cut plate from the corresponding duct. See [Figure 35](#).

**Figure 35.** Duct for the Coupling with a TSQ Quantum Mass Spectrometer

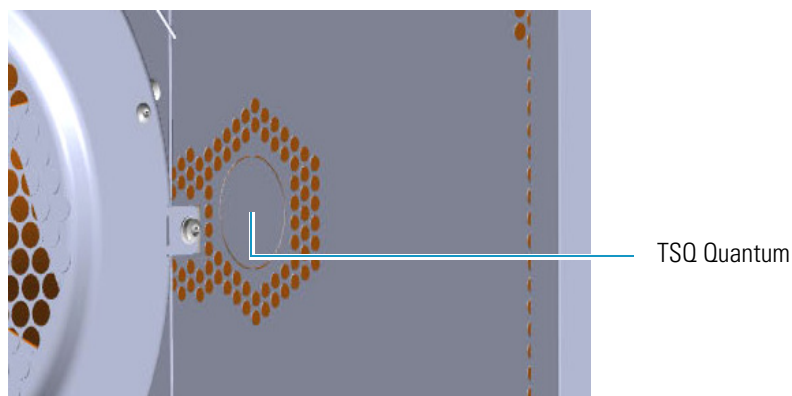


- d. Place the GC right panel for your MS and attach the screw holding it in place.
3. Introduce the transfer line inner tube into the oven through the duct provided.
4. Attach the transfer line to the GC column using the proper nut and ferrule.



**CAUTION - INSTRUMENT DAMAGE:** Make sure that the GC column has been conditioned before connecting it to the transfer line. The material released from the column, (column bleed), during conditioning may contaminate the detector.

5. Close the front door of the GC.
6. Connect the transfer line heater to the **Aux Temperature/Cryo Module**. See [Figure 36](#).

**Figure 36.** Transfer Line MS Heaters

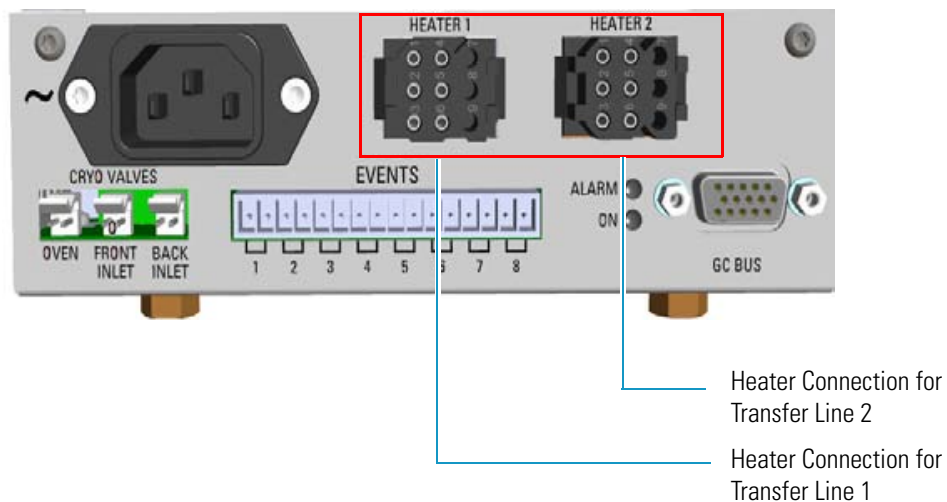


7. To tune and set the TSQ Quantum working conditions, refer to the relevant *User Guide* and *Hardware Manual*.

## Making the Duct to Couple with a DFS, IRMS, or ICP-MS Mass Spectrometer

### ❖ To make the duct for the transfer line


1. Remove the left/right side panel according to the high resolution mass spectrometer to couple.
  - a. Open the front door of the GC.
  - b. Use a T20 Torxhead screwdriver to loosen the left/right side panel screw from the interior front panel. Save the screw because it will be reused later.
  - c. Slide the panel towards the back of the instrument up to the stop.
  - d. Remove the panel by pulling it outward. Be aware that the ground wire is attached to the panel.
2. Prepare the proper duct for the transfer line inner tube.
  - a. Use a T20 Torxhead screwdriver to loosen the duct plate screws from the left/right panel. See [Figure 37](#).

**Figure 37.** Left/Right Side Panel

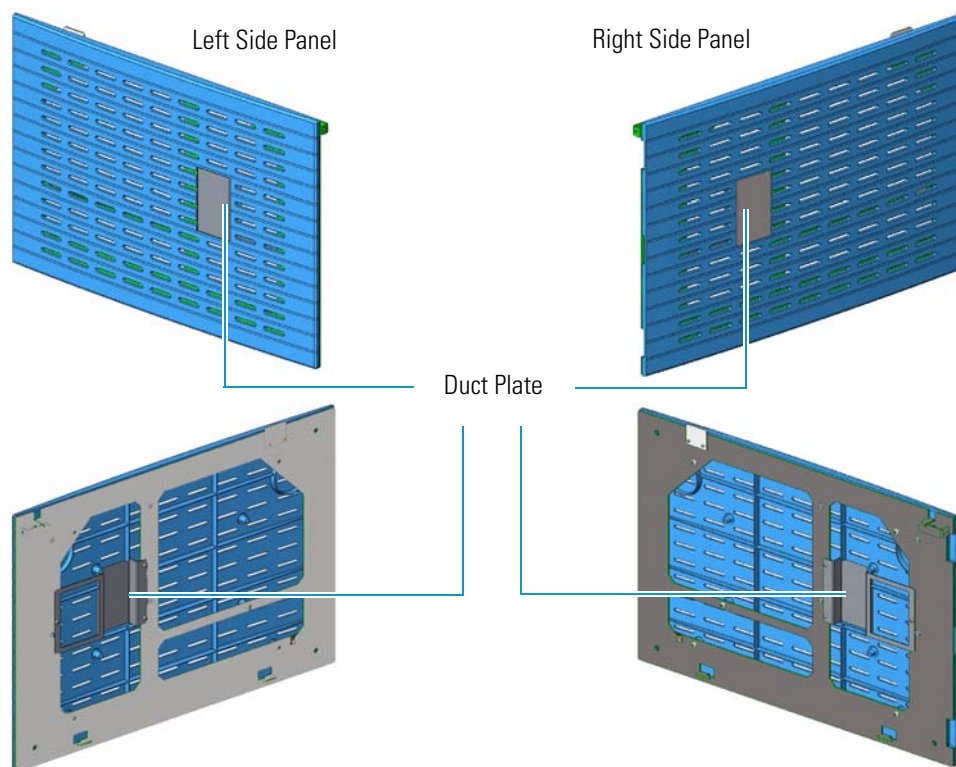

- b. According to the configuration of your high resolution mass spectrometer, on the left/right exterior wall of the oven box, remove the partial cut plate of the duct of interest, and then the pre-shaped plug of insulating material from the duct provided. See [Figure 38](#).


**Figure 38.** Left/Right Exterior Oven Wall Box View

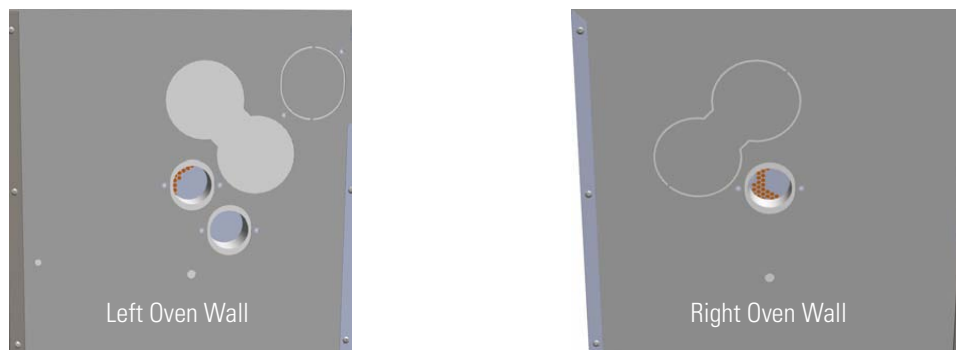

**IMPORTANT** Save the pre-shaped plug of insulating material in a safe place because it could be reused.

- c. On the interior of the oven box remove the partial cut plate from the corresponding duct. See [Figure 39](#).

**Figure 39.** Left/Right Interior Oven Wall Box View

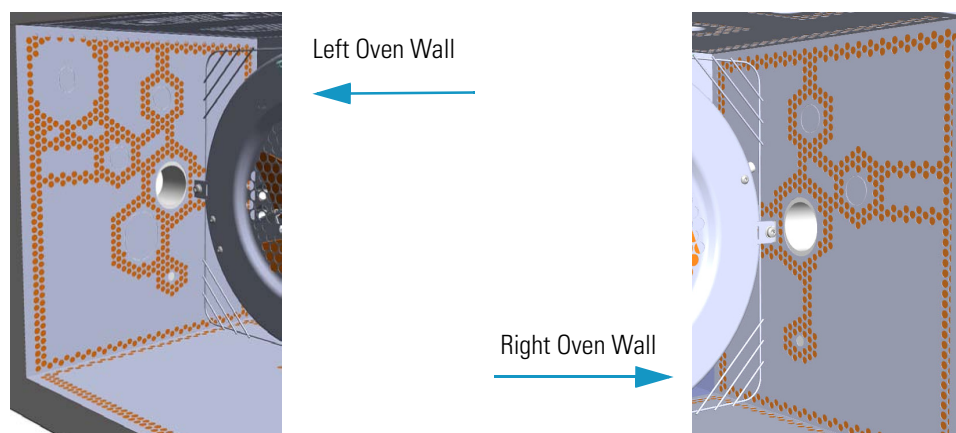


- d. Place the left/right panel for MS and attach the screw holding it in place.
3. Introduce the transfer line inner tube into the oven through the duct provided.
4. Attach the transfer line to the GC column using the proper nut and ferrule.



**CAUTION - INSTRUMENT DAMAGE:** Make sure that the GC column has been conditioned before connecting it to the transfer line. The material released from the column (column bleed) during conditioning may contaminate the detector.

5. Close the front door of the GC.
6. Connect the transfer line heater if required.
  - a. Connect the transfer line heater to the **Aux Temperature/Cryo Module**. See [Figure 40](#).

**Figure 40.** Transfer Line MS Heaters

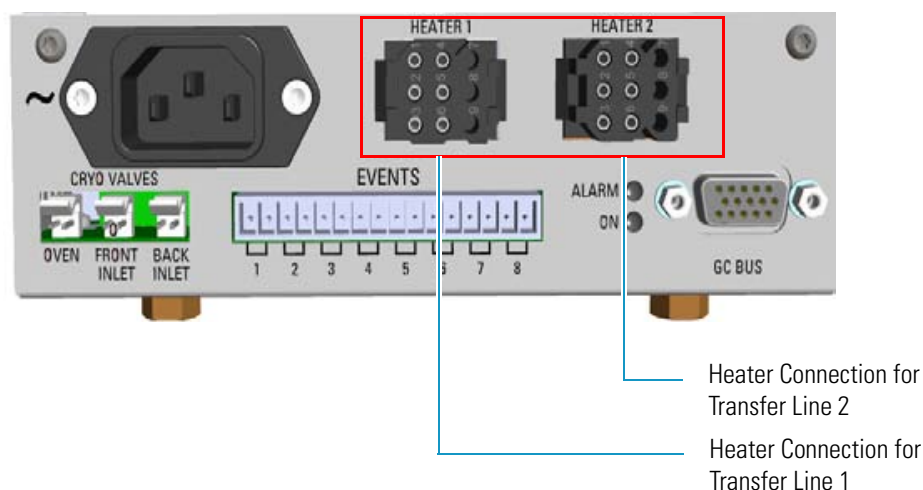

7. Tune and set the DFS, IRMS, and ICP-MS working conditions of the high resolution mass spectrometer, referring to the relevant manuals.

## Setting Handshake Parameters

When your GC is connected to a mass spectrometer, set the handshaking parameters as shown in [Table 1](#).

**Table 1.** GC Handshaking Parameters When A Mass Spectrometer is Connected

| Parameter        | Set To:     |
|------------------|-------------|
| Remote Start In  | High to Low |
| Inhibit Ready In | When High   |
| End of Run Out   | High to Low |
| Start of Run Out | High to Low |
| GC Ready Out     | When Low    |
| Prep Run Out     | When Low    |

### ❖ To configure handshaking parameters through the touch screen

1. In the main menu select the **Configuration** icon, the Configuration menu appears.
2. In the Configuration menu, select the **Handshake** icon to open the relevant submenu.
3. Set the **Handshaking parameters** as required, then return to main menu.

### ❖ To configure handshaking parameters through the Chromatography Data System

1. Launch the Data System. In the relevant Configuration page specify the handshaking parameters.

## Installing the Autosampler

This section provides instruction for installing and connecting a TriPlus RSH, a TriPlus 100 Liquid Sampler, a TriPlus, an AI/AS 1310 - AI/AS 3000 II, or two AI/AS 1310 - AI/AS 3000 II in Gemini configuration.

### Related Topics


- [“Mounting an Autosampler on the GC”](#) on [page 32](#)
- [“Connecting the Autosampler”](#) on [page 34](#)

## Mounting an Autosampler on the GC


To mount the autosampler use the appropriate support and fix it to the installation holes provided on the top of the GC. See [Figure 41](#) and [Figure 42](#).

For further details please refer to the manuals of your sampling system.

**Figure 41.** Installation Holes and Sampler Supports for TriPlus RSH and TriPlus 100 Liquid Sampler

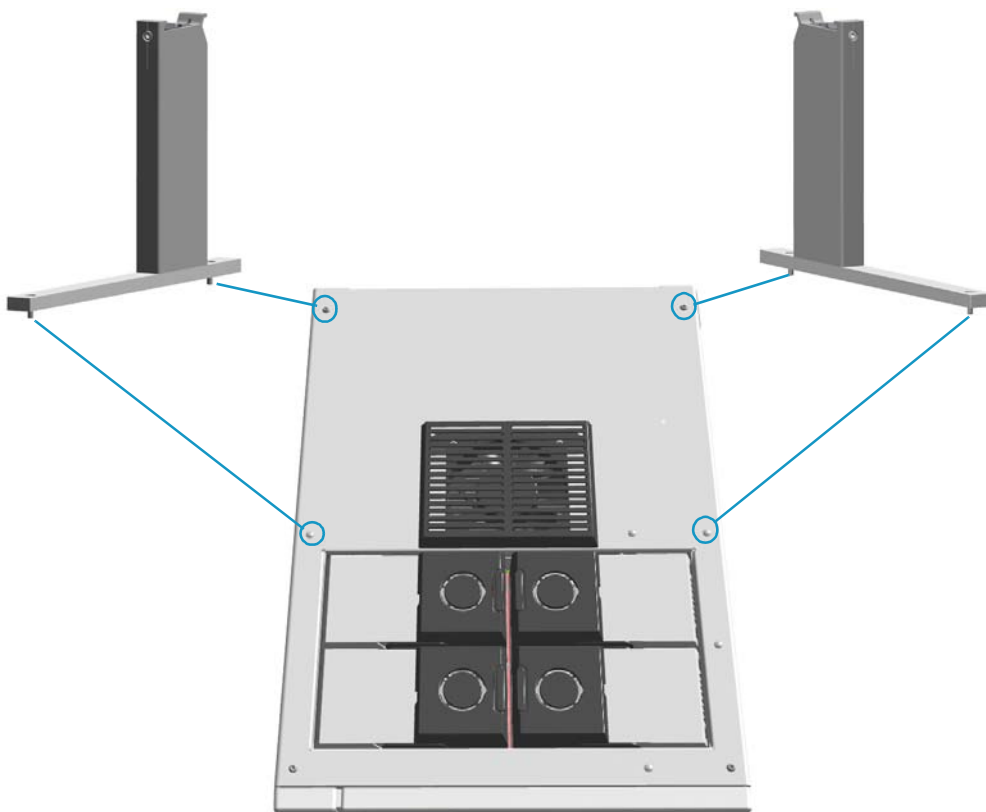

**Figure 42.** Installation Holes and Sampler Supports for AI/AS 1310 - AI/AS 3000 II Autosampler

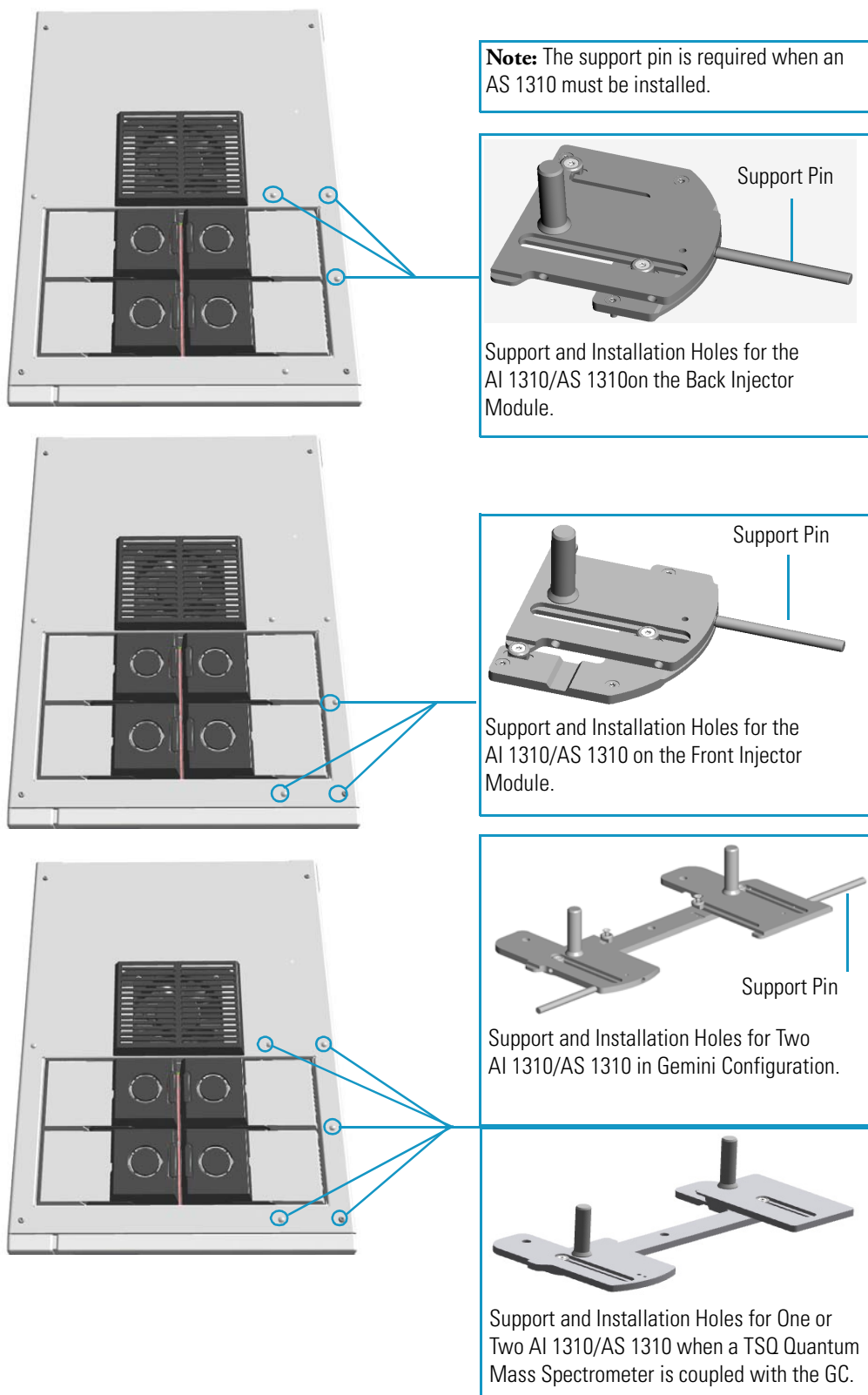




## Connecting the Autosampler

There are two possible cabling patterns for connecting the TRACE 1300/TRACE 1310 to a single autosampler:


- The GC and a TriPlus, AI/AS 1310, or AI/AS 3000 II autosampler are controlled via a single serial cable connected. In this case, a RS-232 connection is also needed between the GC and the autosampler.
- In the Chromatography Data System (CDS) **Serial Port** configuration window, this option is termed **Through GC**. Through GC is the default connection. See [Figure 43](#).

**Figure 43.** Through GC Connection

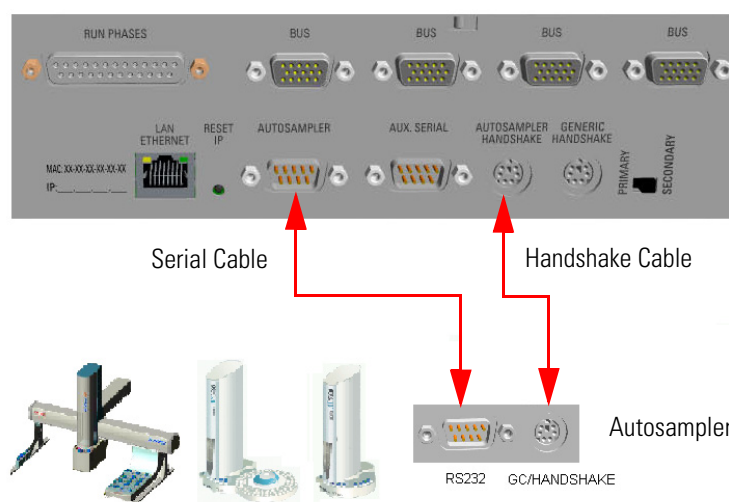

- The GC is connected to the Chromatography Data System (CDS) via LAN cable. The autosampler is connected to the GC via RS-232 serial line and handshake cables. See [Figure 44](#).


**Figure 44.** Separated Connections

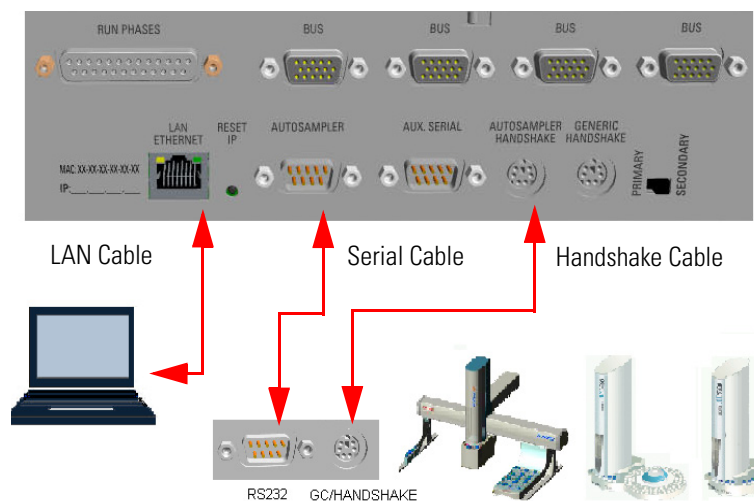




When two AI/AS 1310 or AI/AS 3000 II autosamplers are in use in **Gemini** configuration, perform the hardware setup (**HWSetup 1** or **HWSetup 2**) according to the Thermo Scientific Chromatography Data System in use, and to the High Throughput, Single, or Confirmation operating mode. See [Figure 45](#), [Figure 46](#), and [Table 2](#).

**Figure 45.** Connections of Two AI/AS Autosamplers - HWSetup 1

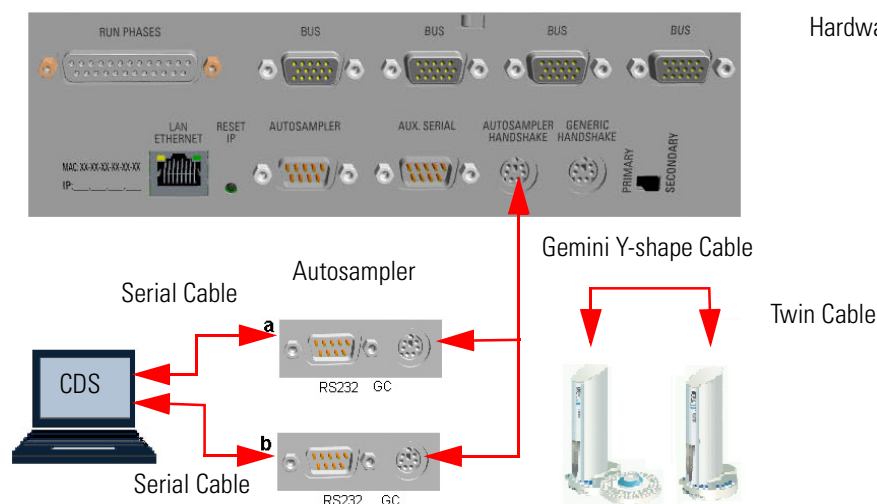

**Figure 46.** Connections of Two AI/AS Autosamplers - HWSetup 2

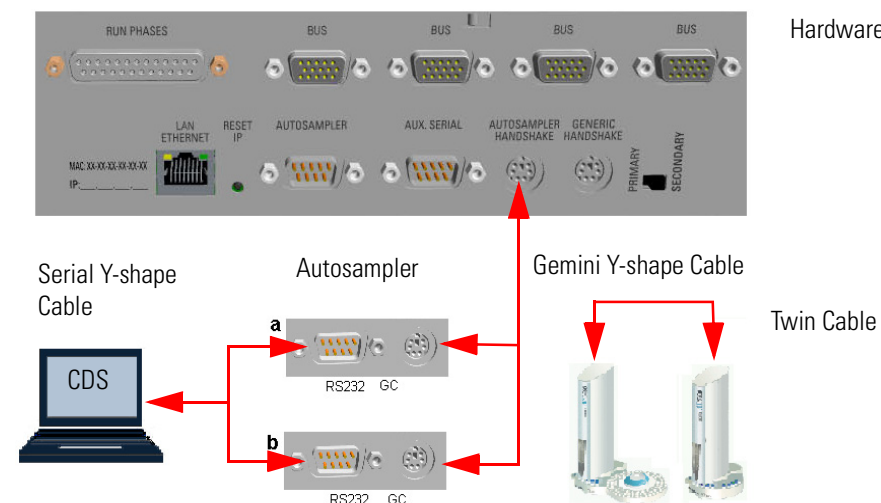

**Table 2.** Hardware Setup Comparison Table (Sheet 1 of 2)

| Chromatography Data System (CDS) | High Throughput Mode                       | Single Mode <sup>1</sup>                                | Confirmation Mode                          |
|----------------------------------|--|---|--|
| Chromeleon 7.2                   | HWSetup 1<br>See <a href="#">Figure 49</a> | HWSetup 1<br>See <a href="#">Figure 49</a>              | HWSetup 2<br>See <a href="#">Figure 49</a> |
| Xcalibur                         | NA <sup>2</sup>                            | HWSetup 2 <sup>1</sup><br>See <a href="#">Figure 50</a> | HWSetup 2<br>See <a href="#">Figure 50</a> |

**Table 2.** Hardware Setup Comparison Table (Sheet 2 of 2)

| Chromatography Data System (CDS) | High Throughput Mode                         | Single Mode <sup>1</sup>                                | Confirmation Mode                          |
|----------------------------------|--|---|--|
| TraceFinder 3.2                  | NA <sup>2</sup>                              | HWSetup 2 <sup>1</sup><br>See <a href="#">Figure 50</a> | HWSetup 2<br>See <a href="#">Figure 50</a> |
| Chrom-Card 2.12                  | Alternative<br>See <a href="#">Figure 50</a> | HWSetup 2 <sup>1</sup><br>See <a href="#">Figure 50</a> | HWSetup 2<br>See <a href="#">Figure 50</a> |
| ChromQuest 5; SP2015 R2          | HWSetup 1<br>See <a href="#">Figure 49</a>   | HWSetup 1<br>See <a href="#">Figure 49</a>              | HWSetup 1<br>See <a href="#">Figure 49</a> |
| EzChrom 3.22; SP2015 R2          | HWSetup 1<br>See <a href="#">Figure 49</a>   | HWSetup 1<br>See <a href="#">Figure 49</a>              | HWSetup 1<br>See <a href="#">Figure 49</a> |


<sup>1</sup> The unplugging of the power cord from the unused AI 1310/AS 1310 autosampler could be required

<sup>2</sup> Not Applicable


- The GC and the TriPlus RSH or TriPlus 100 Liquid Sampler autosampler are connected to the Chromatography Data System (CDS) via LAN cable.

The autosampler is connected to the GC through the Y-shape handshake cable. See [Figure 47](#) and [Figure 48](#).

**Figure 47.** Connection of GC and TriPlus RSH/ TriPlus 100 Liquid Sampler via LAN

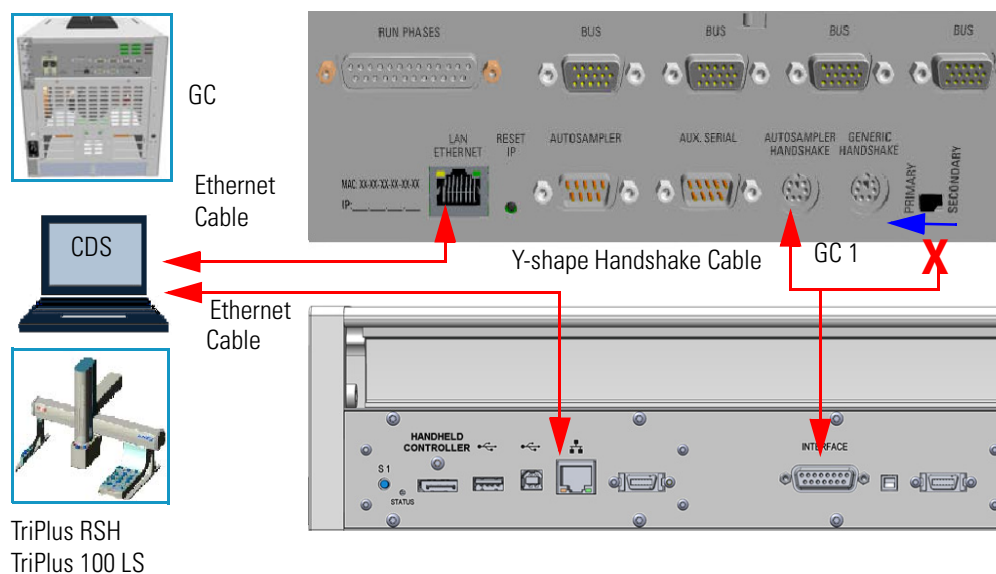

**Figure 48.** Connection of TriPlus RSH/ TriPlus 100 Liquid Sampler to two GC via LAN

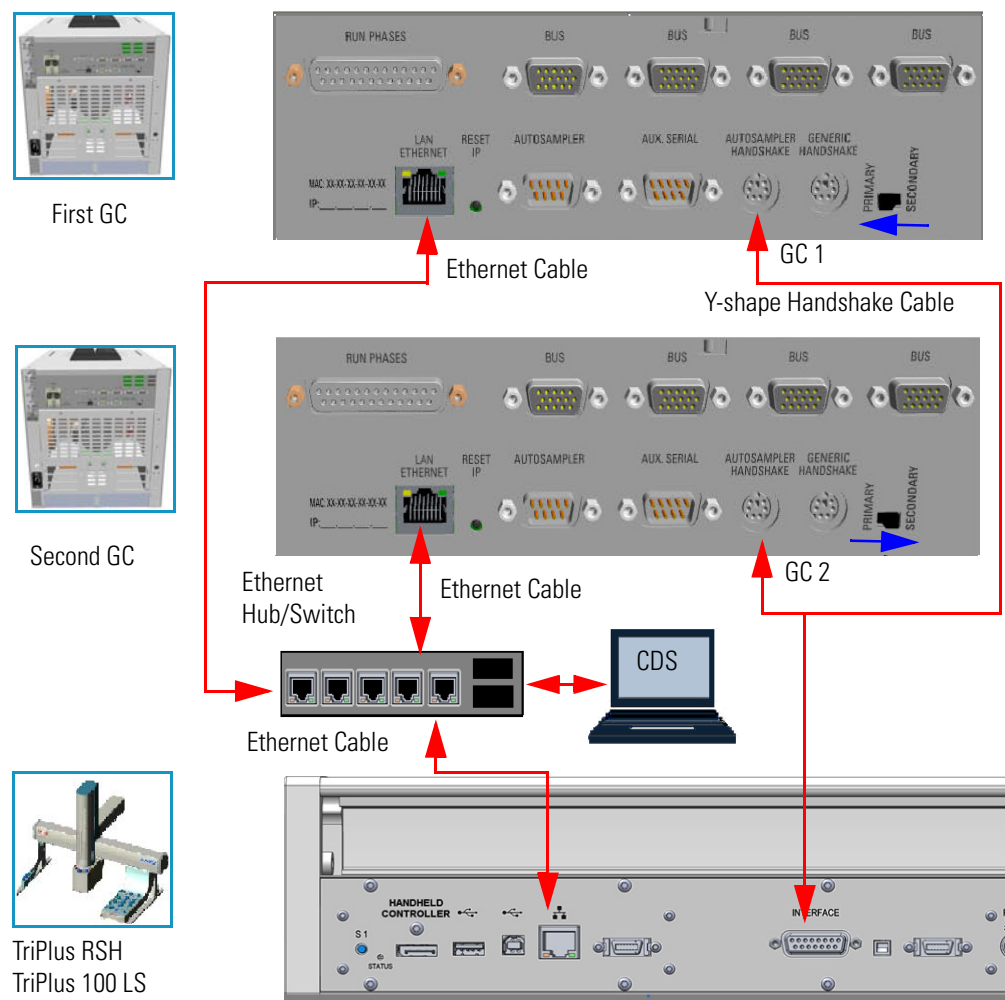

## Connect the Autosampler to the GC System

To connect an autosampler to the GC system see the following procedures:

- “To connect a TriPlus RSH/TriPlus 100 Liquid Sampler autosampler” on page 38
- “To connect a second GC coupled with the same TriPlus RSH/TriPlus 100 Liquid Sampler autosampler” on page 38
- “To connect a a TriPlus autosampler” on page 38
- “To connect a second GC is also coupled with the same TriPlus autosampler” on page 38
- “To connect an AI/AS 1310 - AI/AS 3000 II autosampler” on page 39
- “To perform the HWSetup 1 of two AI/AS 1310 - AI/AS 3000 II in Gemini Configuration” on page 39
- “To perform the HWSetup of two AI/AS 1310 - AI/AS 3000 II in Gemini Configuration” on page 40

❖ **To connect a TriPlus RSH/TriPlus 100 Liquid Sampler autosampler**

1. Use the Ethernet cable to connect the LAN connector on the back of the TriPlus RSH/TriPlus 100 Liquid Sampler Control Board to the computer, and verify that the LAN control LED lights up when connected to a powered-on computer.
2. Using the Y-shape handshake cable provided, connect:
  - a. the main connector to the connection marked **INTERFACE** on the back of the control interface of the TriPlus RSH/TriPlus 100 Liquid Sampler.
  - b. the first branch of the Y-shape handshake cable, labelled **GC 1** to the marked **AUTOSAMPLER HANDSHAKE** on the back panel of the GC.
3. Set the Primary/Secondary switch on **PRIMARY** position.
4. Establish communication between the TriPlus RSH/TriPlus 100 Liquid Sampler and the GC following the instructions detailed in the *TriPlus RSH Hardware Manual* or *TriPlus 100 Liquid Sampler Hardware Manual*.

❖ **To connect a second GC coupled with the same TriPlus RSH/TriPlus 100 Liquid Sampler autosampler**

1. Connect the second branch of the Y-shape handshake cable, labelled **GC 2**, to the connector marked **AUTOSAMPLER HANDSHAKE** on the back panel of the second GC.
2. Set the Primary/Secondary switch on **SECONDARY** position.
3. Establish communication between the TriPlus RSH/TriPlus 100 Liquid Sampler and the GC following the instructions detailed in the *TriPlus RSH Hardware Manual* or *TriPlus 100 Liquid Sampler Hardware Manual*.

❖ **To connect a TriPlus autosampler**

1. Using the cable provided, connect the 6-pin female connector marked **HANDSHAKE MAIN** located on the back portion of the crossrail X and the connector marked **AUTOSAMPLER HANDSHAKE** on the GC back panel.
2. Using the cable provided, connect the 9-pin female connector marked **RS232-1** located on the back portion of the crossrail X to 9-pin connector marked **AUTOSAMPLER 1** on the GC back panel.
3. Connect the power cable to the AC Input connector (Mains socket) on the autosampler and to the wall outlet.

❖ **To connect a second GC is also coupled with the same TriPlus autosampler**

1. Using the cable provided, connect the 6-pin female connector marked **HANDSHAKE VIRTUAL** located on the back portion of the crossrail X to the connector marked **AUTOSAMPLER HANDSHAKE** on the back panel of the second GC.

2. Using the cable provided, connect the 9-pin female connector marked **RS232-2** located on the back portion of the crossrail X to the 9-pin connector marked **AUTOSAMPLER 2** on the GC back panel.
3. Connect the power cable to the AC Input connector (Mains socket) on the autosampler and to the wall outlet.


❖ **To connect an AI/AS 1310 - AI/AS 3000 II autosampler**

1. Using the cable provided, connect the 6-pin female connector marked **GC** located on the sampling unit back side to the connector marked **AUTOSAMPLER HANDSHAKE** on the GC back panel.
2. Using the cable provided, connect the 9-pin female connector marked **RS232** located on the sampling unit back side to the 9-pin connector marked **AUTOSAMPLER 1** on the GC back panel.
3. Connect the power cable to the AC Input connector (Mains socket) on the autosampler and to the wall outlet.

❖ **To perform the HWSetup 1 of two AI/AS 1310 - AI/AS 3000 II in Gemini Configuration**

Figure 49 shows the hardware setup 1 of two AI/AS 1310 or AI/AS 3000 II in Gemini configuration.

**Figure 49.** Gemini Configuration - Hardware Setup 1

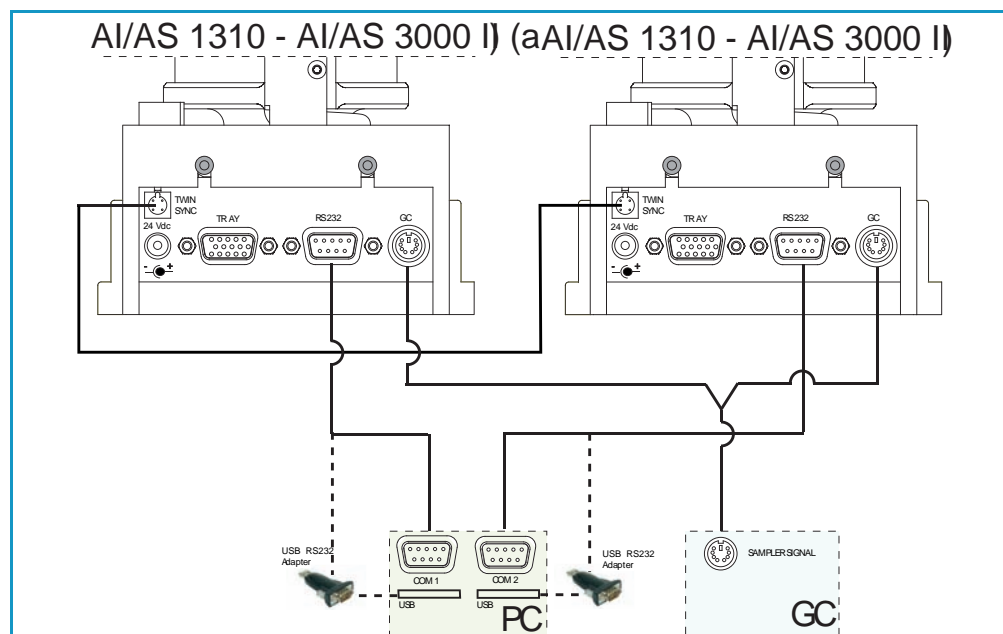


1. Using the cable provided (PN 23043672), connect the 4-pin male connector marked **TWIN SYNC**, on the back of the sampling unit of the sampler “a”, to the 4-pin male connector marked **TWIN SYNC** on the back of the sampling unit of the sampler “b”.
2. Using the **Gemini Y** shape cable provided (PN 23043623), connect:

- a. the 6-pin main connector labelled **GC** to the connector marked **AUTOSAMPLER HANDSHAKE** on the back of the GC.
  - b. the first branch of the cable, labelled **AI/AS 1**, to the 6-pin female connector marked **GC** on the back of the sampling unit of the sampler “a”.
  - c. the second branch of the cable, labelled **AI/AS 2**, to the 6-pin female connector marked **GC** on the back of the sampling unit of the sampler “b”.
3. Using the RS232 cable (PN 23043453) provided in the standard outfit of each AI/AS 1310, connect:
    - a. the 9-pin connector marked **RS232** of the first sampler to a PC serial port (e.g. COM 1), or to an USB port, interposing an USB-RS232 adapter not provided.
    - b. the 9-pin connector marked **RS232** of the second sampler to a second PC serial port (e.g. COM 2), or to an second USB port, interposing a second USB-RS232 adapter not provided.
  4. Connect the power cable to the AC Input connector (Mains socket) on the autosampler and to the wall outlet.

❖ **To perform the HWSetup of two AI/AS 1310 - AI/AS 3000 II in Gemini Configuration**

Figure 50 shows the hardware setup 2 of two AI/AS 1310 or AI/AS 3000 II in Gemini configuration.

**Figure 50.** Gemini Configuration - Alternative Hardware Connection (AHC)

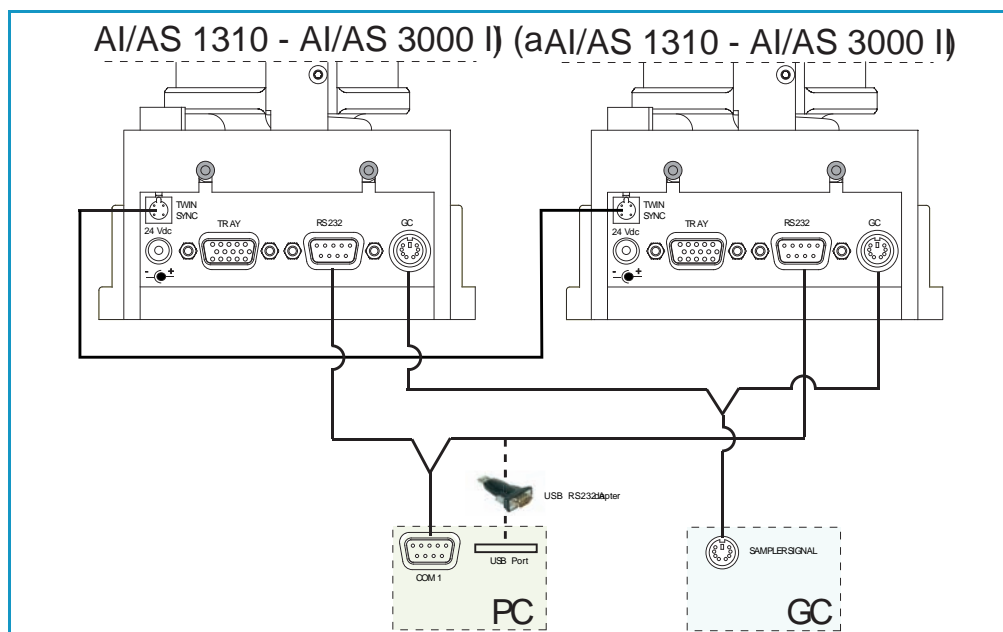

1. Using the cable provided (PN 23043672), connect the 4-pin male connector marked **TWIN SYNC**, on the back of the sampling unit of the sampler “a”, to the 4-pin male connector marked **TWIN SYNC** on the back of the sampling unit of the sampler “b”.

2. Using the **Gemini Y** shape cable provided (PN 23043623), connect:
  - a. the 6-pin main connector labelled **GC** to the connector marked **AUTOSAMPLER HANDSHAKE** on the back of the GC.
  - b. the first branch of the cable, labelled **AI/AS 1**, to the 6-pin female connector marked **GC** on the back of the sampling unit of the sampler “**a**”.
  - c. the second branch of the cable, labelled **AI/AS 2**, to the 6-pin female connector marked **GC** on the back of the sampling unit of the sampler “**b**”.
3. Using the **Serial Y** shape cable provided (PN 23043675), connect:
  - a. the 9-pin main connector to a PC serial port (e.g. COM 1) or to an USB port, interposing an USB-RS232 adapter not provided.
  - b. the first branch of the cable, labelled **AI/AS Primary**, to the 9-pin female connector marked **RS232** on the back of the sampling unit of the sampler “**a**”.
  - c. the second branch of the cable, labelled **AI/AS Secondary**, to the 6-pin female connector marked **RS232** on the back of the sampling unit of the sampler “**b**”.
4. Connect the power cable to the AC Input connector (Mains socket) on the autosampler and to the wall outlet.

## Setting Autosampler Handshake Parameters

When your GC is connected to a AI/AS 1310, AI/AS 3000, TriPlus, TriPlus RSH, or TriPlus 100 Liquid Sampler set the handshaking parameters as shown in [Table 3](#).

**Table 3.** Handshaking Parameters When An Autosampler is Connected

| Parameter        | Set To:     |
|------------------|-------------|
| Remote Start In  | High to Low |
| Inhibit Ready In | Neither     |
| End of Run Out   | High to Low |
| Start of Run Out | High to Low |
| GC Ready Out     | When Low    |
| Prep Run Out     | When Low    |

### ❖ To configure handshaking parameters through the touch screen

1. In the main menu select the **Configuration** icon, the Configuration menu appears.
2. In the Configuration menu, select the **Handshake** icon to open the relevant submenu.
3. Set the **Handshaking parameters** as required, then return to main menu.

❖ **To configure handshaking parameters through the Chromatography Data System**

1. Launch the Data System. In the relevant Configuration page specify the handshaking parameters.

## Installing the Data System Software

The data systems as Chromeleon, Xcalibur, Chrom-Card, or ChromQuest, are designed to be compatible with commercially available computers and requiring the use of Windows™ XP, Windows™ Vista, or Windows™ 7 operating system.

**IMPORTANT** Before installing the Data System, please make sure that any previous version is removed from the PC hard drive.

❖ **To install the data system software**

1. Remove the current version.
  - a. Select **Control Panel > Add/Remove Programs**.
  - b. In the dialog window displayed, select the current Data System software version to remove.
  - c. Click **Add/Remove**.
2. Install the new version.
  - a. Insert the CD/DVD provided and start **Setup.exe**. Click **Next** several times.
  - b. Continue the installation by following the instructions displayed.
  - c. At the end of the installation, reboot the computer if required.
  - d. Start the Data System you installed selecting it in the page **Start-Program** of the desktop.
  - e. The **Main Page** of the data system will be displayed.


For details about the data system in use, please refer to the relevant *User Guide*.



## Making Power Connections

### INSTRUMENT DAMAGE

A label on the electronic module indicates the 120 Vac or 230 Vac power supply. If your power supply line does not match the power supply required, DO NOT CONNECT AND DO NOT POWER ON THE GC.

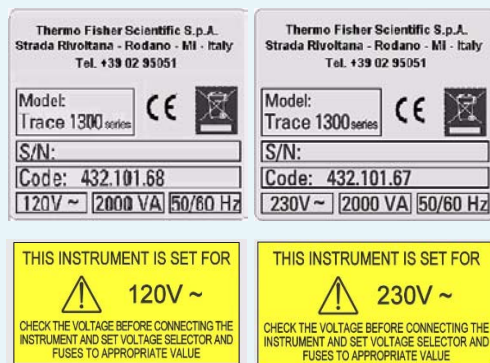


**IMPORTANT** The external modules must be powered at the same line voltage of the main GC system.

### ❖ To connect the power to the GC

1. Connect the power cable to the AC Input connector (Main socket) on the GC, and to the wall outlet. See [Figure 51](#).

**Figure 51.** GC Power Connection

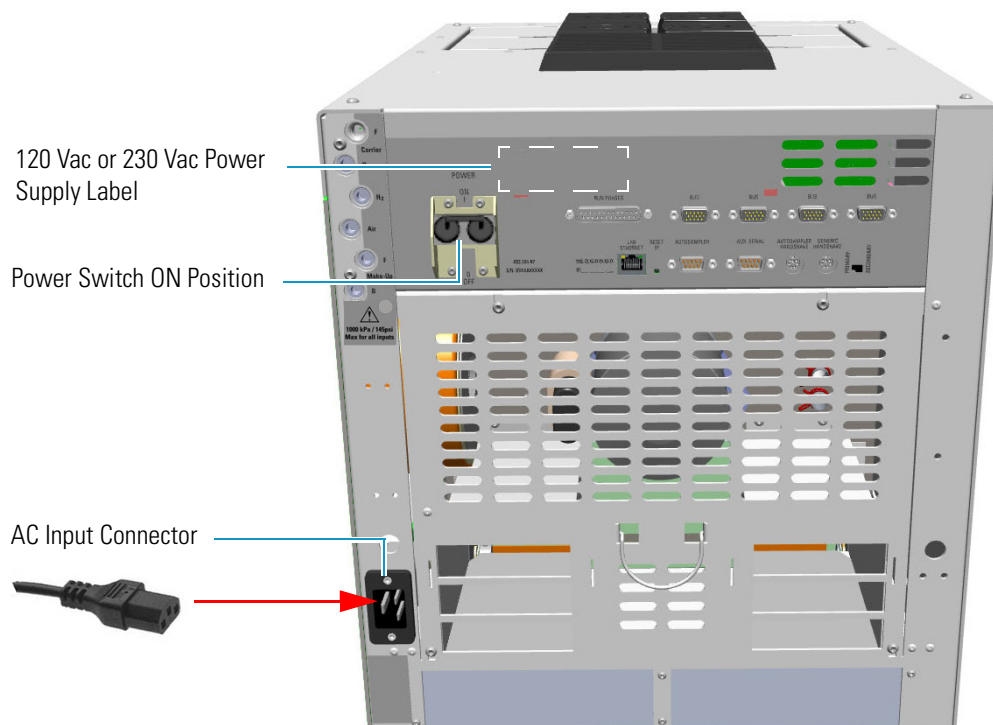


2. If external modules as **Aux Temperature/Cryo Modules**, **NPD Thermionic Source Power Module**, or both are present, connect the power cable to the AC Input connector on the front of each external module, and to the wall outlet. See [Figure 52](#).



**IMPORTANT** The external modules must be powered at the same line voltage of the main GC system.

**Figure 52.** External Module Power Connection

Aux Temperature/Cryo Module

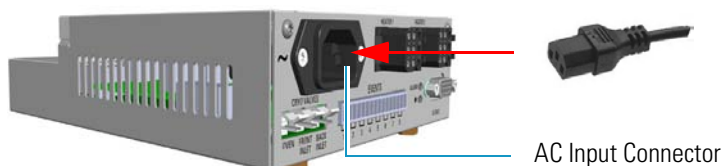

NPD Thermionic Source Power Module



- **Aux Temperature/Cryo Module** — The LED marked **On** lights up solid green after the GC power on.
- **NPD Thermionic Source Power Module** — The LED marked **On** blinks green at the plug-in of the power cable, and becomes solid green at the power on of the thermionic source.

**Note** The modules **Auxiliary Gas**, **Generic Detector Interface**, and **Analog Outputs Interface**, receive the power supply from the electronic module of the GC through the **GC Bus** connection.



The LED marked **On** lights up solid green after the GC power on.

## Setting the LAN Communication


This section provides instructions for setting the desired IP and the LAN communication port for the TRACE 1300/TRACE 1310. See “[Making the LAN Setup](#)” on [page 46](#). Before making the LAN setup, please read the following note.

**Note** The TRACE 1300/TRACE 1310 allows LAN (Local Area Network) capabilities by the presence of an RJ45 connector on the Electronic Module, and two LEDs for LAN activity.

- **IP Address** — The GC is shipped with a factory IP address, which may not match the LAN needs of the sites where the instrument will be installed. To change the default values, contact your LAN administrator and ask for the IP address to be assigned, the netmask, and eventually the port.
  - The IP address is a 3 digits x 4 fields number given by the network administrator - e.g. 192.168.127.10
  - The netmask is a 3 digits x 4 fields number given by the network administrator - e.g. 255.255.255.0
  - The port is a number given by the network administrator - for example 2551
- **Reset Button** — To reset the LAN parameters (IP address, communication port, and so on) to the default values, insert a small screwdriver into the hole, then push the reset button for at least five seconds.
- **Network Cables** — Two network cables are included in the standard outfit:
  - a reversed **RED** patch for the computer to TRACE 1300/TRACE 1310 direct connection. This is used for the initial setup operation and when the instrument is not connected in remote mode.
  - a standard **GREY** patch for the TRACE 1300/TRACE 1310 to local area network connection.

[Figure 53](#) shows the LAN schematic connection between the TRACE 1300/TRACE 1310 and the PC.

**Figure 53.** LAN Schematic Connections

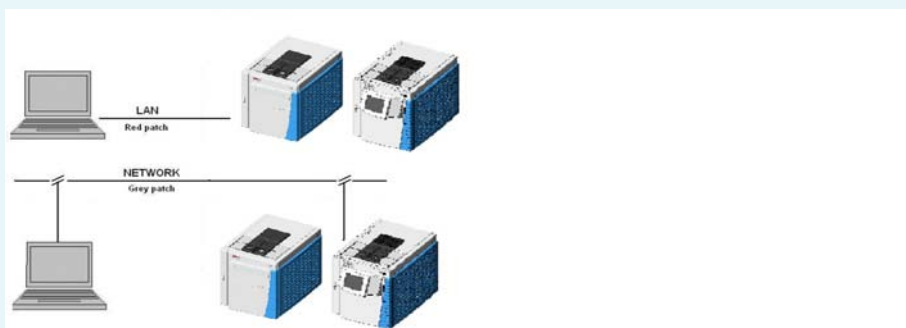


## Making the LAN Setup

### ❖ To set the LAN communication

1. Using the LAN RED patch, connect a desktop or laptop PC directly to the RJ45 connector marked **LAN/Ethernet** on the Electronic Module of the GC. See [Figure 54](#).

**Figure 54.** GC/PC LAN Connection

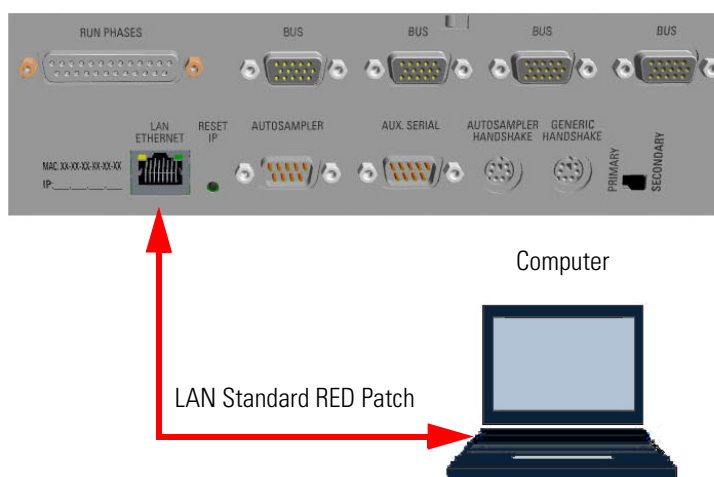

2. Power on the GC by placing the power switch in the On (up) position marked **I**.
3. Power on the PC.
4. Start the **TS1300ST** program to begin the installation on the PC. Follow the instructions step by step to complete the installation.

**Note** The **TS1300ST** program is used for setting the **LAN** Communication parameters, and for updating the **Firmware** version on the GC.




Updating the **Firmware** version on the GC must be carried out by a Thermo Fisher Scientific authorized technical personnel.

Use <http://www.gc-gcms-customersupport.com/WebPage/Share/Default.aspx> address to contact your local Thermo Fisher Scientific office or affiliate GC-GC/MS Customer Support


5. Run the **TS1300ST** program; the initial page appears and you will be prompted to enter the password. See [Figure 55](#).

**Figure 55.** Initial Page and Password

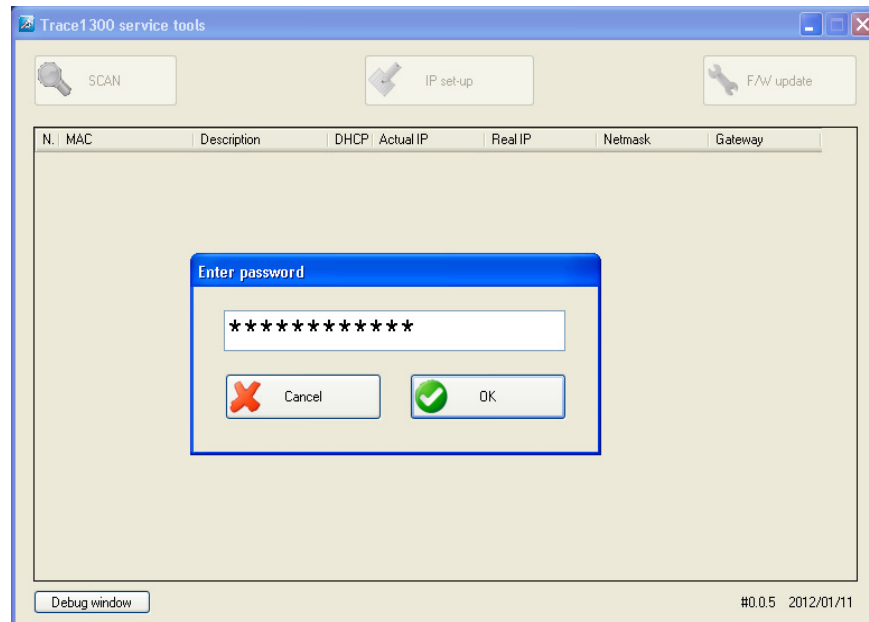


6. In the text box enter the password **ThermoFisher**, then press **OK**.
7. The system starts scanning the equipment powered on and connected to the LAN. During this stage, the three buttons on the top of the page are disabled. See [Figure 56](#).

**Figure 56.** Research of the Equipment

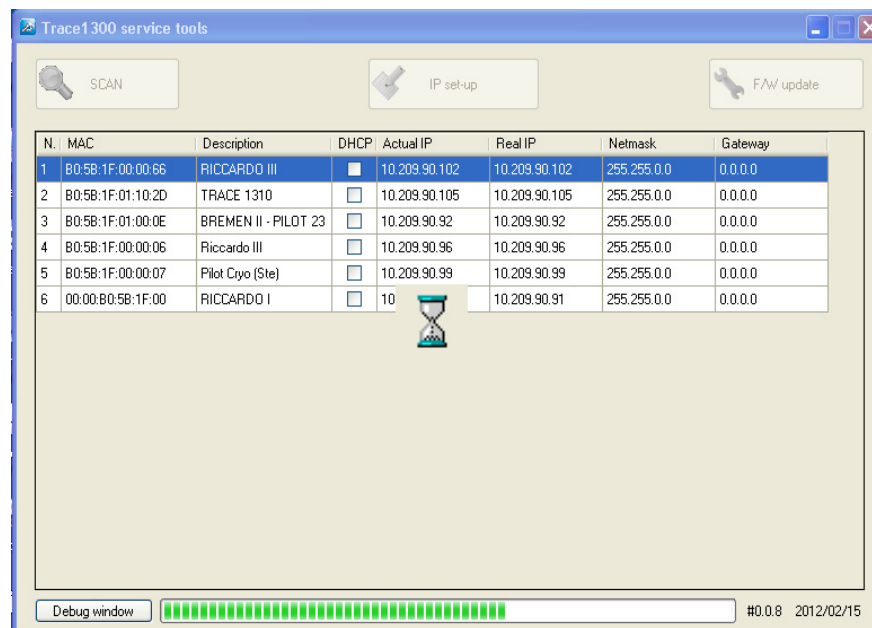

8. At the end of the scan, all the equipment powered on and connected to the LAN is displayed. Only the **SCAN** button is enabled. See the example in [Figure 57](#).


Figure 57. Equipment Connected to the LAN

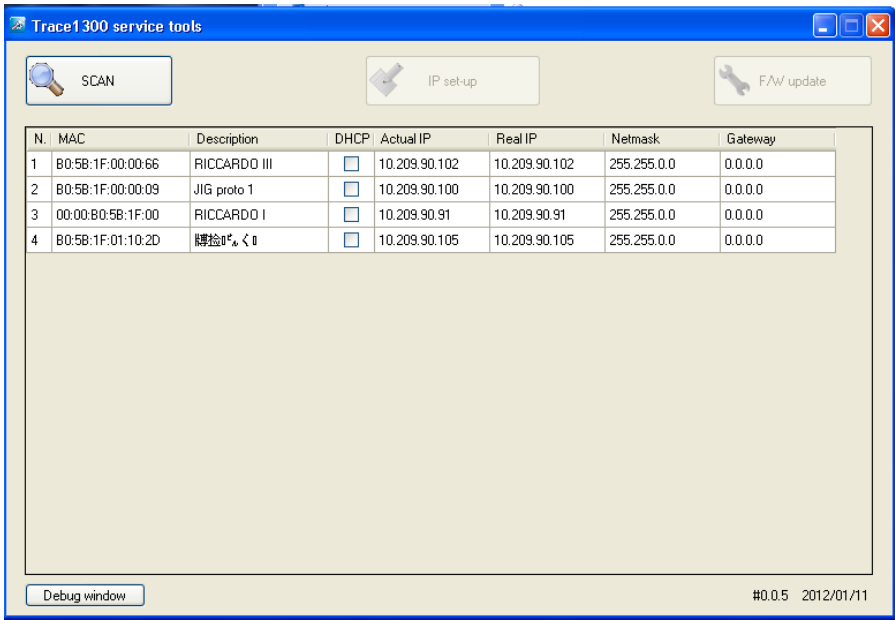

**Note** To repeat the scan, click **SCAN**.

9. Select the GC of interest to set up. All the three buttons on the top of the page are enabled. See the example in [Figure 58](#).

Figure 58. Selection of the GC of Interest

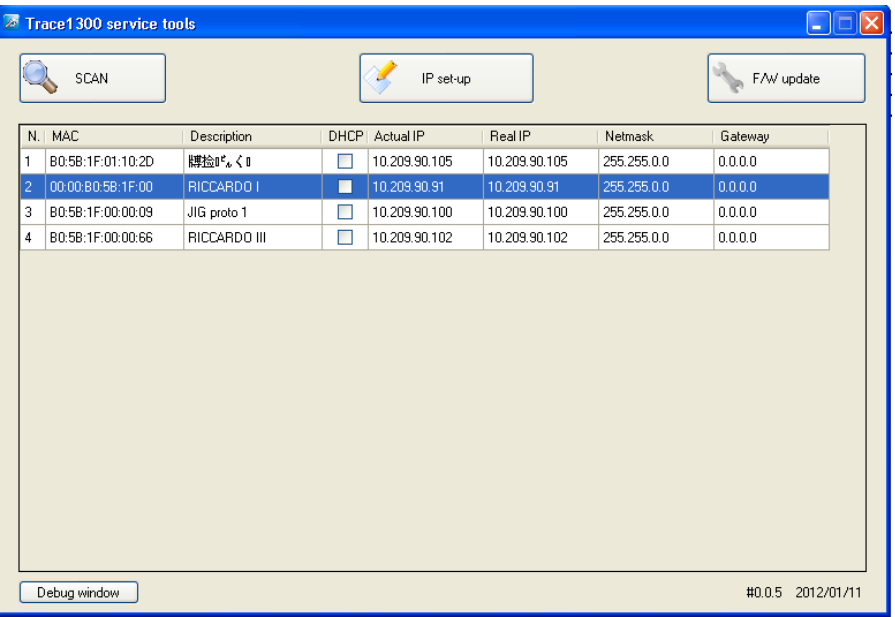

10. Click **IP set-up**. The **Remote Settings** page is displayed. See the example in [Figure 59](#).

**Figure 59.** Remote Settings Page

The RemoteSettings dialog box contains the following fields and controls:


- MAC: 00:00:B0:5B:1F:00
- Description: RICCARDO I
- DHCP: ☐
- IP: 10.209.90.91
- Netmask: 255.255.0.0
- Gateway: 0.0.0.0
- Buttons: Cancel (with a red X icon) and OK (with a green checkmark icon).

The editable parameters are: **Description**, **DHCP**, **IP**, **Netmask**, and **Gateway**.

After changing parameters, press **OK** to confirm the new settings.

11. Close the TS1300ST program by pressing  on the right upper corner of the page, or selecting the command **Close** in the drop-down menu. See [Figure 60](#).

**Figure 60.** Close Commands

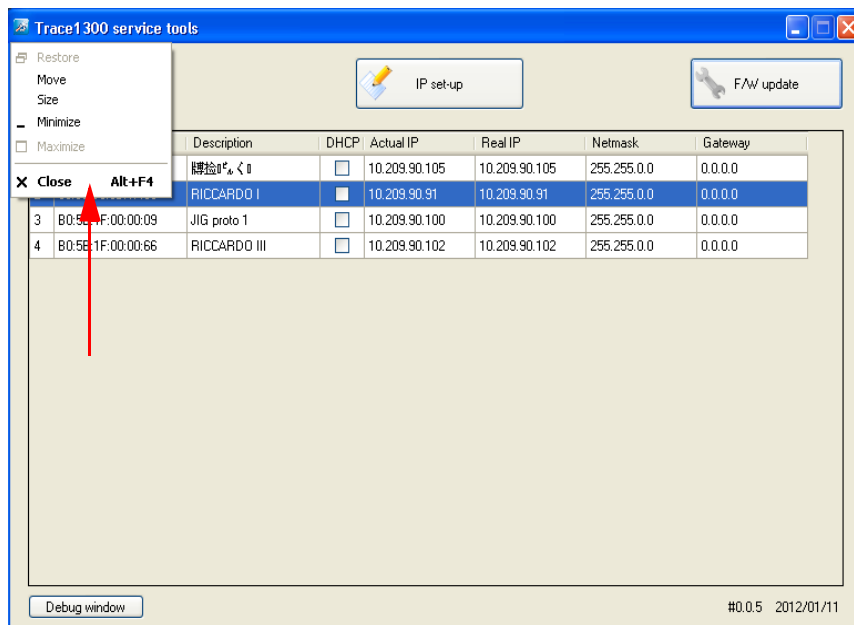

12. Configure the Data System (Xcalibur, Chrom-Card, ChromQuest, or Chromeleon) in use. See the section [“Configuring the Data System”](#) on [page 50](#).

## Configuring the Data System

You should use the self adhesive labels found in the standard outfit of the Data System to annotate the IP address and the TCP Port that are set in the instrument. Stick the label where it can be easily found when necessary to configure the data system.

1. Launch the Data System. In the relevant Configuration page specify a direct TCP/IP address.
2. Enter the IP address of the GC and set the socket used.

Standard operation parameters are set by default; however, your GC may be hidden behind a Firewall that may prevent Port 2551 from being used.

You can set an alternative port number. The number of the socket entered in this box must correspond to the port assigned to the instrument setup.



**IMPORTANT** If the connection is performed through hubs over a 10 Mbit/s network, it is suggested that no more than five GC units be connected on the same network trunk. If you have a switched network, this warning does not apply.

## Column Installation Requirements

The column is where the separation takes place. It should be positioned inside the oven on the column rack. The column ends should align correctly with the injector and the detector bottom fittings. The GC oven controller accurately controls the column temperature. Each column has a maximum recommended operating temperature. To protect the column from excessively high temperatures, remember to set the **Maximum temperature** parameter for the oven.

**Note** For detailed information on maximum operating temperature, refer to the *column manufacturer's instructions* provided with the column.

### Related Topics

- “Using the Correct Fittings” on page 50
- “Installing the Adapters for Encapsulated Graphite Ferrules” on page 52
- “Installing the Column Rack” on page 56


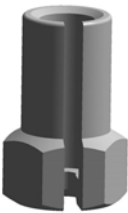
## Using the Correct Fittings

To connect a capillary column to the injector and detector base body, you must use the proper column ferrules and retaining nuts.

- **Retaining Nut** — Dedicated retaining nuts are required to connect capillary columns to injector and detector base bodies, and the **Tee** connector. See [Table 4](#).


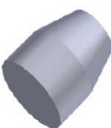



**Table 4.** Column Retaining Nuts for Injectors and Detectors

| Retaining Nut Type   | Injector/Tee Connector              | Detector                |
|--|-------------------------------------|-------------------------|
|  1/4" Hexagonal retaining nut | SSL, SSLBKF                         | FID, NPD, TCD, ECD, FPD |
|  M4 Split retaining nut       | PTV, PTVBKF, SSLBKF (Tee Connector) |                         |

- **Ferrules** — Graphite ferrules, graphite Vespel® ferrules, and encapsulated graphite ferrules are used for many column connections. See [Table 5](#).

**Table 5.** Ferrules for Injectors and Detectors

| Ferrule Type  | Injector/Tee Connector              | Detector                | Temperature Limit |
|---|-------------------------------------|-------------------------|-------------------|
|  Graphite ferrule  | SSL, SSLBKF                         | FID, NPD, TCD, ECD, FPD | 450 °C            |
|  Graphite Vespel® ferrule  | SSL, SSLBKF                         | FID, NPD, TCD, ECD, FPD | 350 °C            |
|  Encapsulated graphite ferrule (Graphpack)<br>See the <b>Note</b> below. | PTV, PTVBKF, SSLBKF (Tee Connector) |                         | 450 °C            |

**Graphite ferrules** are a soft material that is porous to oxygen, making them suitable for most applications except GC/MS interface connections. These easy-to-use ferrules form a soft grip with the column and provide a stable seal. Graphite ferrule are re-usable.

**Graphite/Vespel® ferrules** have a long lifetime and are compatible with GC/MS. These ferrules form a strong grip with the column and cannot be reused as they form a permanent seal with the column. They have a temperature limit of 350 °C, but must be re-tightened after initial temperature cycles.

**Note** An optional adapter kit must be installed if you want to connect the column to the injector and detector by using the encapsulated graphite ferrule. See [“Installing the Adapters for Encapsulated Graphite Ferrules”](#) on page 52.

Table 6, Table 7, and Table 8 list the size of the ferrule to use depending on the internal diameter of precolumn and capillary column.

**Table 6.** Graphite Ferrules Size ID

| Ferrules Size ID | Required for:                      |
|------------------|------------------------------------|
| 0.4 mm           | 0.1 - 0.32 mm ID capillary column  |
| 0.8 mm           | 0.45 - 0.53 mm ID capillary column |

**Table 7.** Graphite Vespel Ferrules Size ID

| Ferrules Size ID | Required for:                       |
|------------------|-------------------------------------|
| 0.4 mm           | 0.1/0.2/0.25 mm ID capillary column |
| 0.5 mm           | 0.32 mm ID capillary column         |
| 0.8 mm           | 0.53 mm ID capillary column         |

**Table 8.** Encapsulated Graphite Ferrules Size ID

| Ferrules Size ID | Required for:                    |
|------------------|----------------------------------|
| 0.4 mm           | 0.25 mm ID capillary column      |
| 0.5 mm           | 0.32 mm ID capillary column      |
| 0.8 mm           | 0.45/0.53 mm ID capillary column |



**CAUTION** Using the wrong size ferrules causes leaks and contamination.




When connecting columns, consider that overtightening compression ferrules does not necessarily produce a stronger, leak-free joint. In fact, very often the reverse is true. Too much pressure can cause a leak in the joint and make it very difficult to reseal the joint when changing columns.

## Installing the Adapters for Encapsulated Graphite Ferrules

This section provides instructions for installing the optional adapters for connecting the capillary column to the SSL injector and FID, NPD, TCD, ECD, and FPD detector by using the encapsulated graphite ferrules (Graph-pack).

The adapters are provided with the kit PN 19050759. The kit includes the parts listed in [Table 9](#).

**Table 9.** Adapters and 1/2-in. Wrench
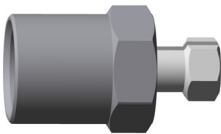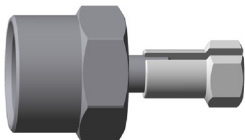
| Part  | Description   |
|---|---|
|  | Adapter for connection to the SSL injector                        |
|  | Adapter for connection to the FID, NPD, ECD, TCD, or FPD detector |
|  | 1/2-in. wrench  |


❖ **To install the adapters for encapsulated graphite ferrules**

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Remove the analytical column.
  - a. Open the front door of the GC.
  - b. Loosen the retaining nut from the injector and detector fitting on the upper interior wall of the GC oven.
  - c. Remove the analytical column with its nut and ferrule from the bottom of the injector and the detector.
5. Remove the bottom parts of the SSL injector. See [Figure 61](#).

**Figure 61.** SSL Injector Bottom Parts Removal

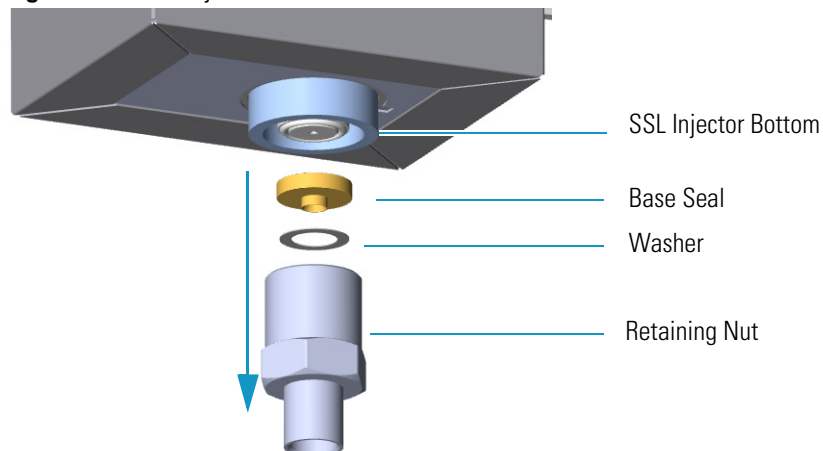

- a. Using the 1/2-in. wrench provided, unscrew the retaining nut with the base seal and the washer from the bottom of the injector.




**CAUTION** Make sure that the liner does not come from the bottom of the injector.

Save the bottom parts of the SSL injector in a safe place because will be reused when you restore the original configuration.


6. Install the adapter on the bottom the injector. See [Figure 62](#).

**Figure 62.** Adapter Installation on SSL Injector Bottom

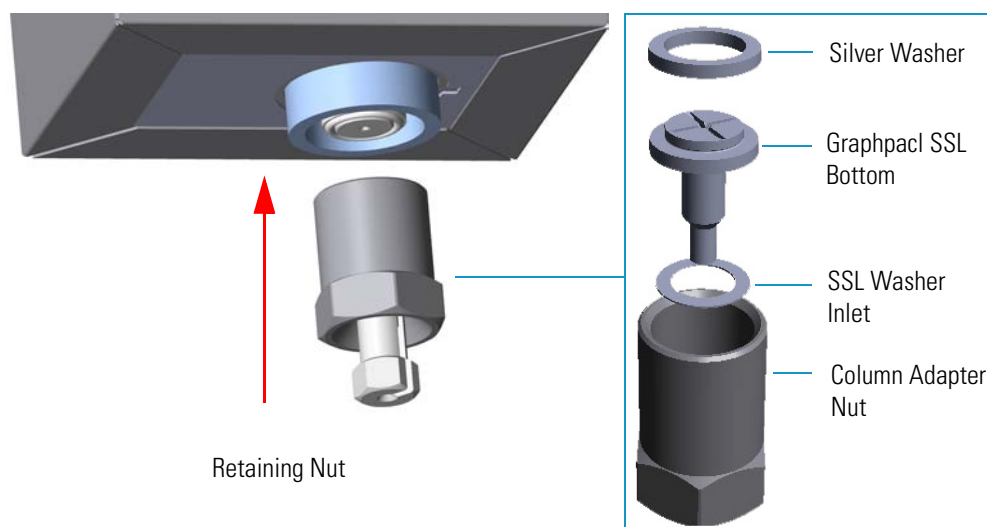

- a. Fingertighten the retaining nut of the adapter until it start to grip the bottom of the injector.
  - b. Use the 1/2-in. wrench to firmly tighten the retaining nut of the adapter with the base seal and the washer from the bottom of the injector.
7. Remove the bottom parts of the FID, NPD, ECD, TCD, or FPD detector. See [Figure 63](#).


**Figure 63.** Detector Bottom Parts Removal

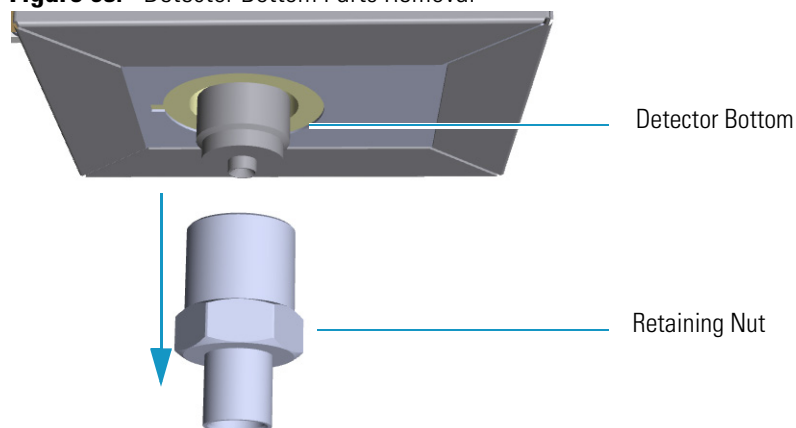

**Note** Save the bottom parts of the detector in a safe place because they will be reused when you restore the original configuration.


- a. Using the 1/2-in. wrench provided, unscrew the retaining nut from the bottom of the detector.
8. Install the adapter on the bottom the detector. See [Figure 64](#).

**Figure 64.** Adapter Installation on the Bottom of the Detector

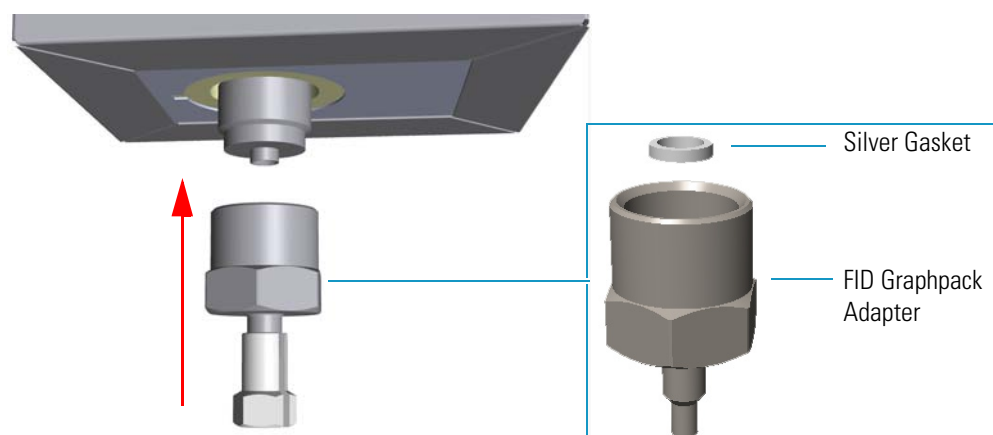

- a. Fingertighten the retaining nut of the adapter until it starts to grip the bottom of the detector.
  - b. Use the 1/2-in. wrench to firmly tighten the retaining nut of the adapter with the base seal and the washer from the bottom of the injector.
9. Connect the column.
- a. Before starting, use a 1/4-5/16-in. wrench to remove the column nuts and ferrules from the base of the injector and detector adapters. See [Figure 65](#).

**Figure 65.** Column Nut and Ferrules

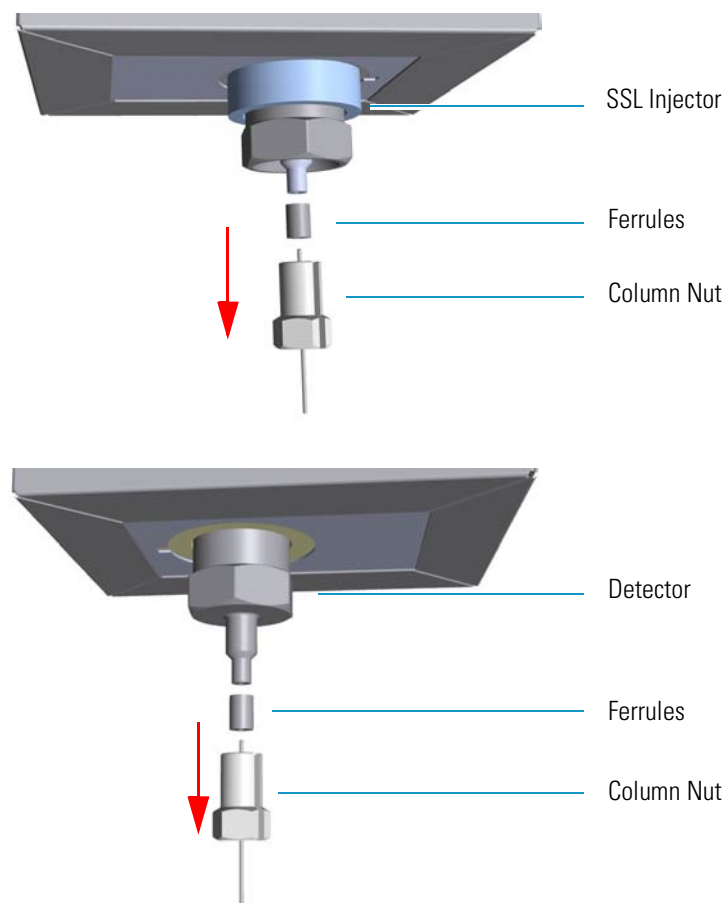


- b. Continue the installation of the capillary column following the instructions in the section [“Installing the Column the First Time”](#) on [page 57](#).

## Installing the Column Rack



**WARNING** Before starting the installation, make sure that the GC is powered off and the power cable is disconnected.

Before connecting the capillary column, make sure that the column rack is installed in the oven. See [Figure 66](#).

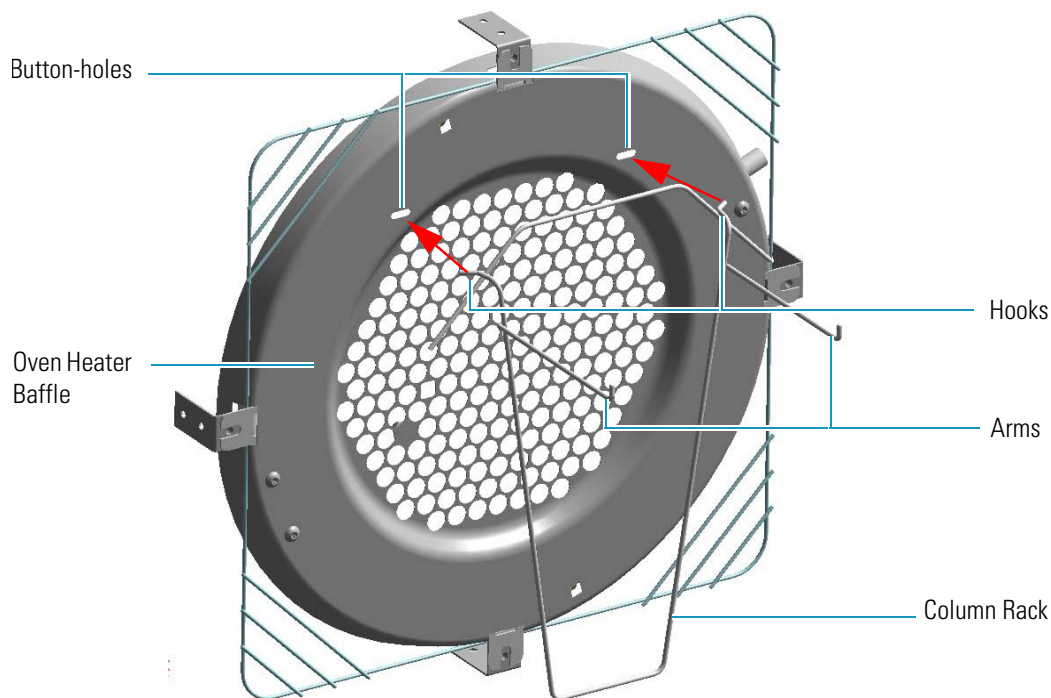**Figure 66.** Column Rack

❖ **To install a column rack**

1. Open the front door of the GC.
2. Slightly press the sides of the column rack and insert the two hooks into the corresponding button-holes on the oven heater baffle.

## Installing the Column the First Time

This section contains instruction for installing the capillary column into the GC oven, for connecting the column ends to the injector and the detector and for carrying out the **column conditioning**, the **leak test** and the **column evaluation**.

**CAUTION** The injector and detector modules are shipped with a plug.




Before installing the column, remove the plug from the injector and detector modules.

Unscrew the column retaining nut from the bottom of the injector and the detector modules, then remove the blind ferrule.

❖ **To install a new column the first time**

1. Install the column into the GC oven.
  - a. Open the front door of the GC.
  - b. Place the column on the two arms of the rack. See [Figure 67](#).

**Figure 67.** Column Installation

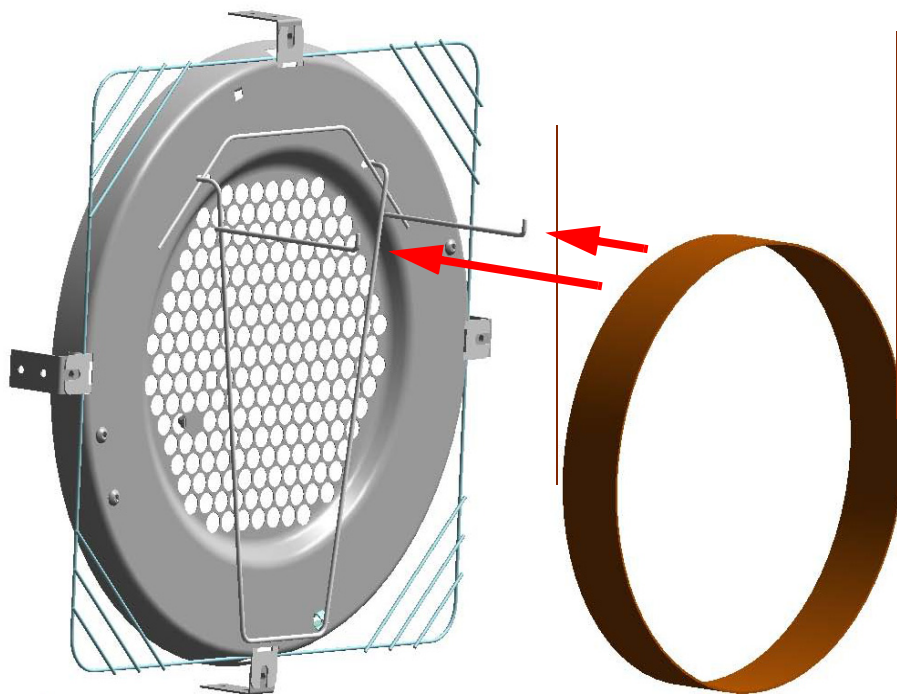

2. Connect the column to the injector inside the oven.

- a. Unwind the column enough to easily connect its ends to the injector and the detector.

**Note** Wear clean, lint- and powder-free gloves when you handle the column and injector ferrule.

- b. Wipe about 100 mm (4 in.) of the column with a tissue soaked in methanol.
- c. Insert the column through the injector retaining nut and the proper ferrule (open end up). Wipe the column again with a tissue soaked in methanol.
- d. Use a scoring wafer to score and break the column about 1 cm (0.4 in.) from the end. Use a magnifying glass to check for an even, flat cut. Repeat if necessary.

**Tip** Slide a notched septum on the column before the injector retaining nut to make it easier to measure the proper distance between the nut and end of the column.

- e. Position the column so that the end of the column extends the proper distance above the end of the ferrule as reported in [Table 10](#).




**Table 10.** Column Insertion Depth For SSL, SSLBKF, HeS-S/SL, PTV, PTVBKF, and GSV Injectors

| Injector         | Column Insertion Depth   |
|------------------|--|
| SSL and HeS-S/SL | <ul style="list-style-type: none"> <li>• 5 mm (splitless)</li> <li>• 5 mm (Helium saver)</li> <li>• 10 mm (split)</li> </ul>                           |
| SSLBKF           | <ul style="list-style-type: none"> <li>• 5 mm (splitless)</li> <li>• 10 mm (split)</li> </ul>  |
| PTV              | <ul style="list-style-type: none"> <li>• 30 mm</li> <li>• As far as possible into the bottom when the PTV is used as an On-Column injector.</li> </ul> |
| PTVBKF           | <ul style="list-style-type: none"> <li>• 30 mm</li> </ul>  |
| GSV              | <ul style="list-style-type: none"> <li>• Insert the column as far as goes and withdrawn about 2 -3 mm</li> </ul>                                       |

- f. Insert the notched septum on the column to hold the retaining nut at this position. Thread the retaining nut into the injector but do not tighten.
  - g. Adjust the column position so that the septum contact the bottom of the retaining nut.
  - h. Finger-tighten the retaining nut until it starts to grip the column plus a quarter turn.
  - i. Remove the notched septum from the column.
3. Open the gas supplies.
  4. Power on the GC.
    - a. Connect the power cable to the AC Input connector on the back of the GC and to the wall outlet.
    - b. Flip up the power switch (breaker) to the position I.
  5. Setup the GC.
    - a. Set the oven and injector temperature to 50 °C.
    - b. Use the column-flowmeter connector to verify that there is flow through the column. If you do not have a flowmeter, dip the column outlet in a small vial of methanol. Bubbles indicate there is flow through the column. If there is no flow, check that the carrier gas is on, the GC inlet is pressurized, and the column is not plugged. If there is still no flow, consult the section **Analytical Troubleshooting** in the *TRACE 1300/TRACE 1310 User Guide*, or contact the Technical Support.
    - c. Allow the column to purge for few minutes.
  6. Perform a column leak check.


- a. Carefully push the capillary column end into the column section of the column-flowmeter connector. See [Figure 68](#).

**Figure 68.** Using a Flowmeter for Leak Check

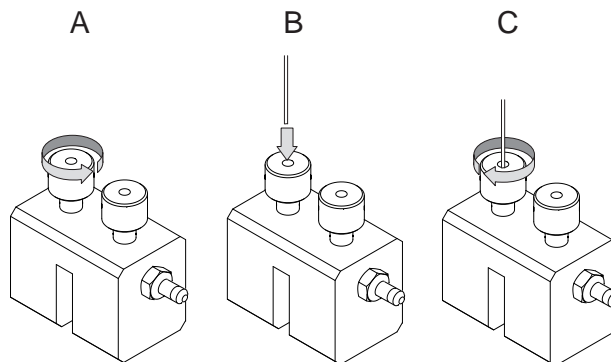

- b. If your GC is equipped with the touch screen as user interface, select the **Leak Check** icon in the **Maintenance** menu, otherwise perform the Leak Check through the Chromatography Data System by selecting the proper function.
  - c. Start the leak check to begin operation. The split and purge valves of the selected channel are automatically closed and the channel is pressurized with carrier gas to the leak check setpoint.
  - d. The system monitors the pressure for one minute. If the pressure does not drop more than the maximum allowed sensitivity value, then the leak check will pass.
  - e. If the leak check did not pass, you should use the leak detector to find and fix the leaks.
- Tip** Leaks can be caused by not tightening the fitting on the column-flowmeter connector. We recommend that you check that fitting before looking elsewhere.
- f. Repeat the leak check until no leaks are indicated.
7. Calibrate the carrier gas flow (column evaluation).
- a. Carefully push the capillary column end into the flow meter section of the column-flowmeter connector. See [Figure 69](#).

**Figure 69.** Column Flowmeter Connector

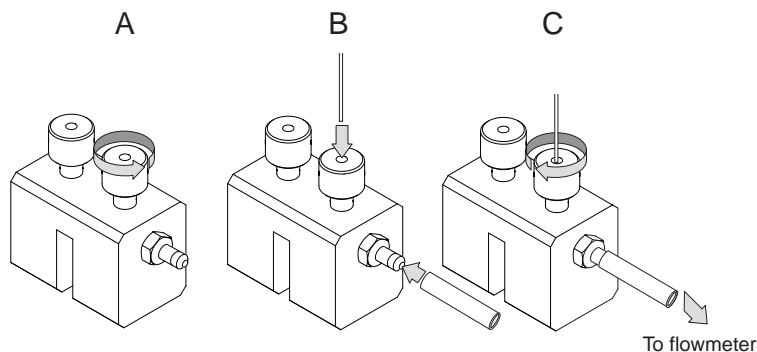

- b. Connect the flowmeter to the dedicated fitting on the column-flowmeter connector.
- c. If your GC is equipped with the touch screen as user interface, select **Back** or **Front Column** icon in the **Configuration** menu, otherwise perform the Column Evaluation through the Chromatography Data System by selecting the proper function.
- d. Select the column and input the physical characteristics of the column.
- e. If a pre-/post-column is present, set the length and nominal internal diameter of the pre-/post-column in the same valid ranges for the column. The following two lines are added to the menu.
- f. According to the physical characteristics of the column, the system calculates and displays the relevant Column K-factor.

**Note** For the most reproducible results, you should conduct a more detailed column evaluation. However, the following steps, while recommended, are not required.

- g. Start column evaluation. At the end of the routine, a message will indicate that evaluation was successful.
  - h. Expect a K-factor of approximately 0.7 – 0.9 for a 15 m, 0.25 mm ID column (1.3 – 2.0 for a 30 m, 0.25 mm ID column). If the column does not report a K-factor within this range or within 0.1 units of the previous stored value, check for a leak or broken column using the leak detector. The K-factor is a measured resistance for the column. A K-factor that is too low may indicate a leak in the system, while a K-factor that is too high may indicate a blockage.
8. Disconnect the column-flowmeter.
    - a. Disconnect the column from the column-flowmeter connector.
    - b. Remove the clear plastic component, including its fittings, from the oven and set it aside.
    - c. Close the GC door.
  9. Condition the column.

The column must be conditioned before inserting it into the detector.

Column conditioning consists of passing a carrier gas flow through the column heated at a programmed temperatures as described in the *column manufacturer's instructions*.

In case the column does not have any column conditioning instructions, perform the column conditioning by setting a final temperature up to 10 °C - 20 °C below its recommended maximum temperature.



**CAUTION** When performing column conditioning, the column should be connected only to the injector leaving the column outlet disconnected to avoid the possibility of contamination of the detector. Do not use hydrogen as the carrier for conditioning! It could vent into the oven and present an explosion hazard.

- a. Run the temperature program that is recommended by the manufacturer.



**INSTRUMENT DAMAGE:** Never exceed the column manufacturer's maximum operating temperature.

10. Connect the column to the detector inside the GC.

- a. Lower the oven temperature to 30 °C and allow it to cool.



**WARNING-BURN HAZARD:** The injector, detector, oven, and transfer line may be hot. Allow them to cool to room temperature before touching them.

- b. Unwind the column enough to easily connect its ends to the injector and the detector..

**Note** Wear clean, lint- and powder-free gloves when you handle the column and injector ferrule.

- c. Wipe about 100 mm (4 in.) of the column with a tissue soaked in methanol.
- d. Use a scoring wafer to score and break the column outlet about 2.5cm (1 in.) from the end. Use a magnifying glass to check for an even, flat cut. Repeat if necessary.
- e. Insert the column through the proper detector retaining nut and ferrule (open end up). Wipe the column again with a tissue soaked in methanol.

**Tip** Slide a notched septum on the column before the detector retaining nut to make it easier to measure the proper distance between the bottom nut and end of the column.

- f. For **FID**, **NPD**, **TCD**, **ECD**, and **FPD**, position the column so that the end of the column extends the proper distance above the end of the ferrule as reported in [Table 11](#). For **PDD** see the instruction described at the [step g](#) on [page 63](#).

**Table 11.** Column Insertion Depth For FID, NPD, TCD, ECD, FPD, and PDD Detectors

| Detector          | Column Insertion Depth                                       |
|-------------------|--|
| FID, NPD, and TCD | Insert the column as far as goes and withdrawn about 2 -3 mm |
| ECD               | 23 mm  |
| FPD               | 125 mm   |
| PDD               | 136 mm   |

- i. For **FID**, **NPD**, and **TCD**, insert the column into the detector, paying attention to not force it further. Finger-tighten the retaining nut, then withdraw the column **2-3 mm**. Tighten the retaining nut an additional a quarter turn.

- ii. For **ECD** and **FPD**, insert the notched septum on the column to hold the retaining nut in this position. Thread the retaining nut into the detector but do not tighten. Adjust the column position so that the septum contact the bottom of the retaining nut.  
Finger-tighten the retaining nut until it starts to grip the column plus a quarter turn.
- iii. Remove the notched septum from the column.
- g. For **PDD** the capillary column must enter **136 mm** into the pre-installed capillary column adapter.
  - i. Make a mark on the column 136 mm from the end.
  - ii. Remove the knurled nut column inlet at the bottom of the detector. Slide the nut over the end of the column, followed by the appropriate column ferrule.
  - iii. Seat the ferrule in the detail of the column adapter and begin sliding the column through the capillary column adapter and into the column inlet.
  - iv. Get the nut started on the threads and tighten it until you feel it contact the ferrule, then back off half a turn.
  - v. Slide the column into the column inlet until the mark is flush with the surface of the knurled nut, and secure the column in the adapter by tightening the knurled nut finger tight only.

**Note** When inserting the capillary column into the PDD detector it might rarely happen to feel a slight resistance. In this case, for proper column installation, pull the column out slightly and adjust the angle before inserting it further.



**IMPORTANT** To install a packed column, the pre-installed capillary column adapter must be replaced with the **packed columns adapter** that enters into the PDD cell for the correct length.

- 11. End of the column installation.
  - a. Close the front door of the GC.



## Performing Routine Maintenance

This chapter provides instructions for performing routine maintenance on the TRACE 1300 and TRACE 1310.

### Contents

- [Read Me First](#)
- [Maintenance Supplies and Tools](#)
- [Maintenance Button](#)
- [Powering On the TRACE 1300/TRACE 1310](#)
- [Shutting Down the TRACE 1300/TRACE 1310](#)
- [Cleaning the Instrument Externally](#)
- [Replacing a Column](#)

# Read Me First

The instrument will be generally serviced by Thermo Fisher Scientific authorized technical personnel for all the warranty period or, after warranty, possibly according to a Programmed Service Contract. For more information contact your local Thermo Fisher Scientific office.

**WARNING** If, for technical reasons, it is necessary to work on parts of the machine that may involve hazardous operations (moving parts, components under voltage, and so on). Thermo Fisher Scientific authorized Technical Support must be called.

This situation can be identified because the access to these moving parts is possible only using a particular tool, and because the concerned removable protective covers bear a warning symbol that draws the operator's attention to the specific warnings included in the documentation accompanying the instrument. In case the work must be carried out by the operator, the latter must prove to be adequately trained to perform the specific maintenance operation.

There are only a few TRACE 1300/TRACE 1310 components that require routine maintenance, depending on the quantity and types of samples you are running. A frequently used instrument will, of course, require more maintenance than an instrument that is rarely used.

- **External Cleaning** — The GC needs to be cleaned when it gets dirty. See the section “[Cleaning the Instrument Externally](#)” on [page 72](#).
- **Column** — You may need to replace the column when your performance degrades and troubleshooting indicates that the column needs maintenance. That may mean that end of the column needs to be trimmed or the column needs to be replaced. See “[Replacing a Column](#)” on [page 73](#).
- **Injector modules** — You may need to install, replace or maintain an injector module. See [Chapter 3, “Performing Injectors Routine Maintenance.”](#)
- **Detector modules** — You may need to install, replace or maintain a detector module. See [Chapter 4, “Performing Detectors Routine Maintenance.”](#)

There are many more components in the TRACE 1300/TRACE 1310 that do not require routine maintenance, but may need to be replaced if there is a problem with the instrument. To replace any component not listed in this chapter, see the following chapters:

- [Chapter 5, “GC Main Frame Advanced Maintenance.”](#)
- [Chapter 6, “Injectors Advanced Maintenance.”](#)
- [Chapter 7, “Detectors Advanced Maintenance.”](#)



**IMPORTANT**

All the maintenance operations must be carried out at low temperature to avoid burns. Therefore, before beginning the sequence, the GC oven, injectors and detectors must be cooled to room temperature, and then the gases supply must be closed.

You can carry out these operations manually, or by pressing the **Maintenance** key. For details refer to the *TRACE 1300 and TRACE 1310 User Guide*.

At the end of the maintenance operations, restore the GC normal working conditions.

**Maintenance of an injector**

Before opening the injector for maintenance, turn the carrier gas off, and wait for the carrier pressure to go to zero.

**When an autosampler is present:**

Move the autosampler away from the injector module to create free space around it:

- If an AI 1310/AS 1310 - AI 3000/AS 3000 II is installed, pull the sampler support plate outwards.
- If a TriPlus RSH, TriPlus 100 Liquid Sampler, or TriPlus is installed, make sure that the standby turret position does not obstacle the injector or detector module to maintain. If yes, change the position of the turret/head.

## Maintenance Supplies and Tools

To perform routine maintenance on the TRACE 1300/TRACE 1310, you will need the following supplies and tools.

- Wrench, open-end, 1/4-in. and 5/16-in., 1/8-in.
- Flathead screwdriver
- 5.5 x 25 mm Slotted Stubby Driver
- Column cutter, wafer (5181-8836, 4/pk)
- T20 Torxhead screwdriver
- T10 Torxhead screwdriver
- T6 Torxhead screwdriver
- 3-mm Allen key wrench
- Electronic flowmeter (Thermo Scientific GFM Pro Flowmeter, or equivalent)
- Electronic leak detector (Thermo Scientific GLD Pro, or equivalent)
- Tweezers (or thin needle-nose pliers) or forceps
- Gloves, heat-resistant (for handling hot parts)

## Cleaning Stainless Steel Components

To clean stainless steel components, you will need:

- Acetone, reagent grade (or other suitable polar solvent)
- Applicators, cotton-tipped
- De-ionized water
- Detergent (Alconox™, Micro®, or equivalent)
- Gas, clean and dry (N<sub>2</sub> or He)
- Gloves, clean, lint- and powder-free, latex or nitrile
- Lint-free cloth
- Ultrasonic cleaner

## Maintenance Button

When pressed, the **Maintenance** button allows cooling at 60 °C of the selected heated zones, which allows the required maintenance. Carrier and detector gases must be closed manually. This function also performs the following actions automatically:


- turns off the flame and the fuel gases of the FID and FPD detectors.
- turns off the thermionic source, the hydrogen, and the air flows of the NPD detector.
- turns off the filaments of the TCD detector.

## Powering On the TRACE 1300/TRACE 1310

### ❖ To power on the GC


1. Install the GC column (see [“Replacing a Column”](#) on [page 73](#)).
2. Open the supply gases.
3. Plug the power cable to the AC Input connector (Main socket) on the GC, and to the wall outlet. See [Figure 70](#).
4. Be sure the carrier gas flowing through the column and the detector gases flowing through the detector.

**Figure 70.** GC Power On

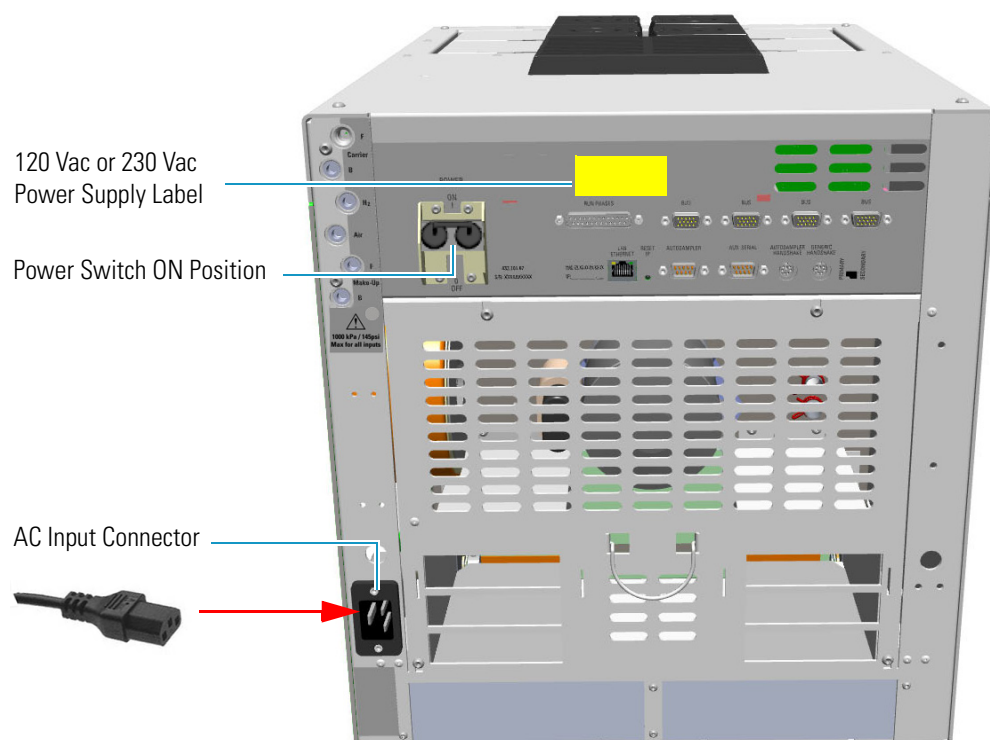

5. If external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.
6. Flip up the power switch, located on the back side of the GC, to the On (up) position marked **I**.
  - a. **TRACE 1300** — When the GC powers on, all the LEDs on the status panel light up simultaneously, afterward the **Power** light becomes a solid green while all the other lights turn off. The GC is now in stand-by status. See [Figure 71](#).

**Figure 71.** TRACE 1300 Status Panel at the GC Power On

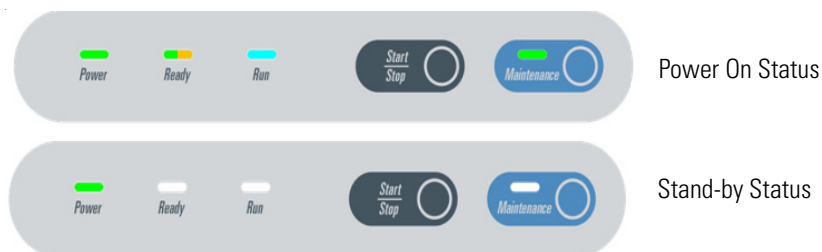

- b. **TRACE 1310** — Check the main menu appears on the touch screen. See [Figure 72](#).

**Figure 72.** TRACE 1310 Touch Screen



7. Open the Thermo Scientific Chromatography Data System installed on the computer.
8. Set the analytical parameters.

## Shutting Down the TRACE 1300/TRACE 1310


### ❖ To shut down the GC

1. Cool down the GC.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

2. If you do not plan to replace the column or perform maintenance on the GC, you do not have to lower the injector temperature.
3. Turn off the carrier gas supply at the tank.
4. Push down the power switch (breaker), located at the back of the instrument, to the position **O** (down). See [Figure 73](#).

**Figure 73.** GC Power Off

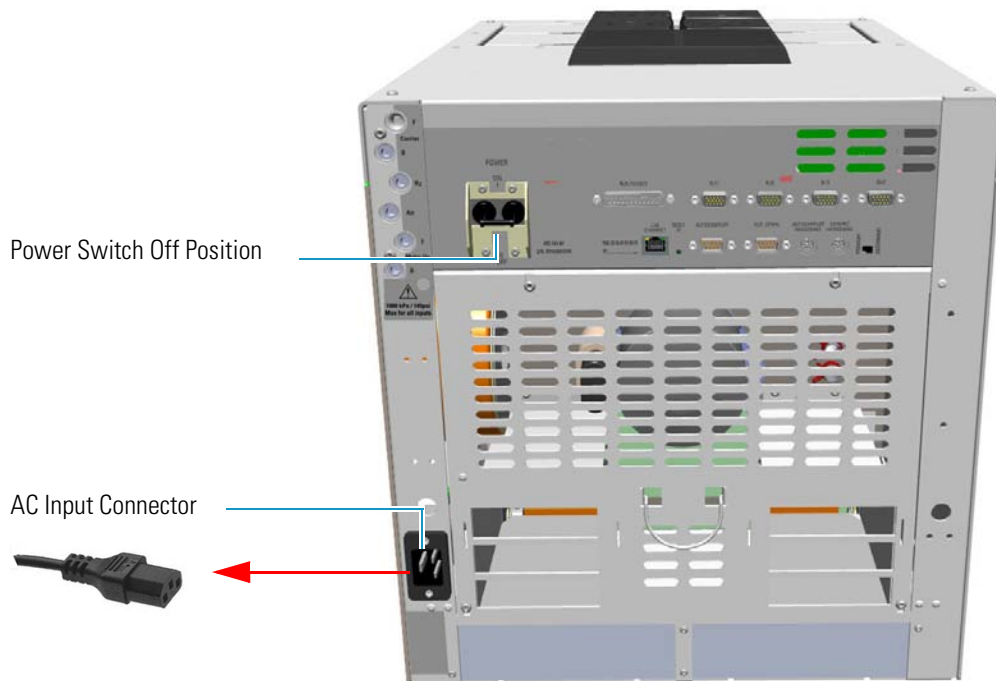

5. If external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.
6. If present, power-off the autosampler by switching off the main power switch, or by unplugging the power cable from the AC input connector, and from the wall outlet.
7. Power-off all the remaining instruments.

## Cleaning the Instrument Externally

Normal usage of the TRACE 1300/TRACE 1310 can cause the exterior to get dirty.

**WARNING** It is your responsibility to avoid that dangerous liquids, materials or both, seeping inside the GC during operation and maintenance.

Solvent must not be used. Do not spray on electrical parts.

### ❖ To clean the instrument externally

1. Place the GC in stand-by condition.
2. Press the **Maintenance** button to cool down the GC, or switch off the heated zones manually.
3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector on the back of the GC, and from the wall outlet.
4. If external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.
5. Externally clean the instrument with a soap and water solution, or with a household non-abrasive product.
  - Pay special attention when cleaning the back side of the instrument. Do not spray, but clean using a cloth imbued with the same substance.
  - Carefully avoid seeping of the products used inside the instrument, particularly when cleaning the grid of the back panel.
  - If you just suspect that a substance used for cleaning or a product submitted to analysis has penetrated inside the instrument, immediately shut down the instrument, and call an authorized customer support engineer for proper actions. The service engineer must be fully informed on the nature of the concerned substance.  
In the event that a hazardous material is spilled on or in the instrument, clean the spill according to the procedures reported in the Material Data Sheet for that substance.
6. Dry with a clean cloth.
7. If external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.
8. Power on the GC.

- a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
- b. Flip up the power switch (breaker) to the position I.
- c. Set the normal injector, detector and GC working conditions.

## Replacing a Column

You might need to replace the column when your performance degrades and troubleshooting indicates that the column needs maintenance. That might mean that the end of the column needs to be trimmed, or the column needs to be replaced.



**WARNING-BURN HAZARD:** The injector and the oven could be hot. Cool to room temperature before touching them.

### ❖ To replace a column

**Note** If you are running samples, stop the acquisition.

1. Remove the current column.
  - a. Press the **Maintenance** button to cool down the GC.
  - b. Set off the carrier and detector gases of the channel of interest, and wait for the carrier pressure to go to zero.
  - c. Open the front door of the GC.
  - d. Unscrew the injector and detector nuts, and remove the column.
  - e. Remove the column from the column rack, and from the GC.
2. Install the new column.
  - a. Place the new column on the two arms of the column rack inside the oven.
3. Connect the new column to the injector inside the GC.
  - a. Unwind the column enough to easily connect its ends to the injector and the detector..

**Note** Wear clean, lint- and powder-free gloves when you handle the column and injector ferrule.

- b. Wipe about 100 mm (4 in.) of the column with a tissue soaked in methanol.
- c. Insert the column through the proper injector retaining nut and ferrule (open end up). If the M4 retaining nut is used, slide it on the column through the side cut. Wipe the column again with a tissue soaked in methanol.
- d. Use a scoring wafer to score and break the column about 1 cm (0.4 in.) from the end. Use a magnifying glass to check for an even, flat cut. Repeat if necessary.

**Tip** Slide a notched septum on the column before the injector retaining nut to make it easier to measure the proper distance between the nut and end of the column.

- e. Position the column so that the end of the column extends the proper distance above the end of the ferrule as reported in [Table 12](#).

**Table 12.** Column Insertion Depth For SSL, SSLBKF, HeS-S/SL, PTV, PTVBKF, and GSV Injectors

| Injector         | Column Insertion Depth   |
|------------------|--|
| SSL and HeS-S/SL | <ul style="list-style-type: none"> <li>• 5 mm (splitless)</li> <li>• 5 mm (He saver)</li> <li>• 10 mm (split)</li> </ul>                               |
| SSLBKF           | <ul style="list-style-type: none"> <li>• 5 mm (splitless)</li> <li>• 10 mm (split)</li> </ul>  |
| PTV              | <ul style="list-style-type: none"> <li>• 30 mm</li> <li>• As far as possible into the bottom when the PTV is used as an On-Column injector.</li> </ul> |
| PTVBKF           | <ul style="list-style-type: none"> <li>• 30 mm</li> </ul>  |
| GSV              | <ul style="list-style-type: none"> <li>• Insert the column as far as goes and withdrawn about 2 -3 mm</li> </ul>                                       |

- f. Insert the notched septum on the column to hold the retaining nut at this position. Thread the retaining nut into the injector but do not tighten.
  - g. Adjust the column position so that the septum contact the bottom of the retaining nut.
  - h. Finger-tighten the retaining nut until it starts to grip the column plus a quarter turn.
  - i. Remove the notched septum from the column.
4. Setup the GC parameters.
    - a. Set the oven and injector temperature to 50 °C.
    - b. Use the column flowmeter connector to verify that there is flow through the column. If you do not have a flowmeter, dip the column outlet in a small vial of methanol. Bubbles indicate there is flow through the column. If there is no flow, check that the carrier gas is on, the GC inlet is pressurized, and the column is not plugged. If there is still no flow, consult the Troubleshooting section or contact Technical Support.
    - c. Allow the column to purge for at least 10 minutes.
  5. Perform a column leak check.




- a. If your GC is equipped with the touch screen as user interface, select the **Leak Check** icon in the **Maintenance** menu, otherwise perform the Leak Check through the Chromatography Data System by selecting the proper function.
- b. Start the leak check to begin operation. The split and purge valves of the selected channel are automatically closed and the channel is pressurized with carrier gas to the leak check setpoint.
- c. The system monitors the pressure for one minute. If the pressure does not drop more than the maximum allowed sensitivity value, then the leak check will pass.

If the leak check did not pass, you should use the leak detector to find and fix the leaks.

**Tip** Leaks can be caused by not tightening the fitting on the column flowmeter connector. We recommend that you check that fitting before looking elsewhere.

- d. Repeat the leak check until no leaks are indicated.
6. Calibrate the carrier gas flow (column evaluation).
- a. Carefully push the capillary column end into the flow meter section of the column flowmeter connector. See [Figure 74](#).

**Figure 74.** Column Flowmeter Connector

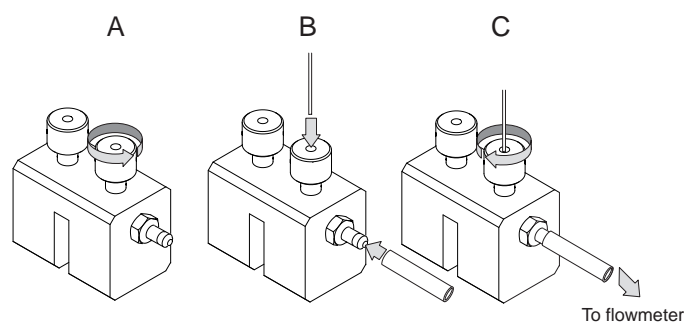

- b. Connect the flowmeter to the dedicated fitting on the column flow meter connector.
- c. If your GC is equipped with the touch screen as user interface, select **Back** or **Front Column** icon in the **Configuration** menu, otherwise perform the Column Evaluation through the Chromatography Data System by selecting the proper function.
- d. Select the column and input the physical characteristics of the column.
- e. If a pre-/post-column is present, set the length and nominal internal diameter of the pre-/post-column in the same valid ranges for the column. The following two lines are added to the menu.
- f. According to the physical characteristics of the column, the system calculates and displays the relevant Column K-factor.

**Note** For the most reproducible results, you should conduct a more detailed column evaluation. However, the following steps, while recommended, are not required.

- g. Start column evaluation. At the end of the routine, a message will indicate that evaluation was successful.
  - h. Expect a K-factor of approximately 0.7 – 0.9 for a 15 m, 0.25 mm ID column (1.3 – 2.0 for a 30 m, 0.25 mm ID column). If the column does not report a K-factor within this range or within 0.1 units of the previous stored value, check for a leak or broken column using the leak detector. The K-factor is a measured resistance for the column. A K-factor that is too low may indicate a leak in the system, while a K-factor that is too high may indicate a blockage.
7. Disconnect the column flowmeter.
- a. Disconnect the column from the column flowmeter connector.
  - b. Remove the clear plastic component, including its fittings, from the oven and set it aside.
  - c. Close the GC door.
8. Condition the column.

The column must be conditioned before inserting it into the detector.

Column conditioning consists of passing a carrier gas flow through the column heated at a programmed temperatures as described in the *column manufacturer's instructions*.

In case the column does not have any column conditioning instructions, perform the column conditioning by setting a final temperature up to 10 °C - 20 °C below its recommended maximum temperature.



**CAUTION** When performing column conditioning, connect the column only to the injector leaving the column outlet disconnected to avoid the possibility of contamination of the detector.

Do not use hydrogen as the carrier for conditioning! It could vent into the oven and present an explosion hazard.

- a. Run the slow temperature program that is recommended by the manufacturer.



**INSTRUMENT DAMAGE:** Never exceed the column manufacturer's maximum operating temperature.

9. Connect the column to the detector inside the GC.
- a. Lower the oven temperature to 30 °C and allow it to cool.



**WARNING-BURN HAZARD:** The injector, detectors, oven, and transfer line may be hot. Allow them to cool to room temperature before touching them.

- b. Unwind the column enough to easily connect its ends to the injector and the detector.

**Note** Wear clean, lint- and powder-free gloves when you handle the column and injector ferrule.

- c. Wipe about 100 mm (4 in.) of the column with a tissue soaked in methanol.
- d. Use a scoring wafer to score and break the column outlet about 2.5cm (1 in.) from the end. Use a magnifying glass to check for an even, flat cut. Repeat if necessary.
- e. Insert the column through the proper detector retaining nut and ferrule (open end up). Wipe the column again with a tissue soaked in methanol.

**Tip** Slide a notched septum on the column before the detector retaining nut to make it easier to measure the proper distance between the bottom nut and end of the column.

- f. For **FID**, **NPD**, **TCD**, **ECD**, and **FPD**, position the column so that the end of the column extends the proper distance above the end of the ferrule as reported in [Table 13](#). For **PDD** see the instruction described at the [step g](#) on [page 77](#).

**Table 13.** Column Insertion Depth For FID, NPD, TCD, ECD, FPD, and PDD Detectors

| Detector          | Column Insertion Depth                                       |
|-------------------|--|
| FID, NPD, and TCD | Insert the column as far as goes and withdrawn about 2 -3 mm |
| ECD               | 23 mm  |
| FPD               | 125 mm   |
| PDD               | 136 mm   |

- i. For **FID**, **NPD**, and **TCD**, insert the column into the detector, paying attention to not force it further. Finger-tighten the retaining nut, then withdraw the column **2-3 mm**. Tighten the retaining nut an additional a quarter turn.
- ii. For **ECD** and **FPD**, insert the notched septum on the column to hold the retaining nut in this position. Thread the retaining nut into the detector but do not tighten. Adjust the column position so that the septum contact the bottom of the retaining nut. Finger-tighten the retaining nut until it starts to grip the column plus a quarter turn.
- iii. Remove the notched septum from the column.
- g. For **PDD** the column must enter **136 mm** into the pre-installed capillary column adapter.
  - i. Make a mark on the column 136 mm from the end.

- ii. Remove the knurled nut column inlet at the bottom of the detector. Slide the nut over the end of the column, followed by the appropriate column ferrule.
- iii. Seat the ferrule in the detail of the column adapter and begin sliding the column through the capillary column adapter and into the column inlet.
- iv. Get the nut started on the threads and tighten it until you feel it contact the ferrule, then back off half a turn.
- v. Slide the column into the column inlet until the mark is flush with the surface of the knurled nut, and secure the column in the adapter by tightening the knurled nut finger tight only.

**Note** When inserting the capillary column into the PDD detector it might rarely happen to feel a slight resistance. In this case, for proper column installation, pull the column out slightly and adjust the angle before inserting it further.



**IMPORTANT** To install a packed column, the pre-installed capillary column adapter must be replaced with the **packed columns adapter** that enters into the PDD cell for the correct length.

10. End of the column installation.
  - a. Close the front door of the GC.

## Performing Injectors Routine Maintenance

This chapter provides instructions for performing routine maintenance on the TRACE 1300 and TRACE 1310 injector modules.

### Contents


- [Maintaining a Split/Splitless Injector \(SSL\)](#)
- [Maintaining a Split/Splitless Injector with Backflush \(SSLBKF\)](#)
- [Maintaining a Gas Sampling Valve Injector \(GSV\)](#)
- [Maintaining an Instant Connect Helium Saver Injector \(HeS-S/SL\)](#)
- [Maintaining a Programmable Temperature Vaporizing Injector \(PTV\)](#)
- [Maintaining a Programmable Temperature Vaporizing Injector with Backflush \(PTVBKF\)](#)

## Maintaining a Split/Splitless Injector (SSL)


This section provides instructions for maintaining a Split/Splitless injector (SSL).

The module and injector components are shown in [Figure 75](#) and [Figure 76](#).

**Figure 75.** SSL Module Components

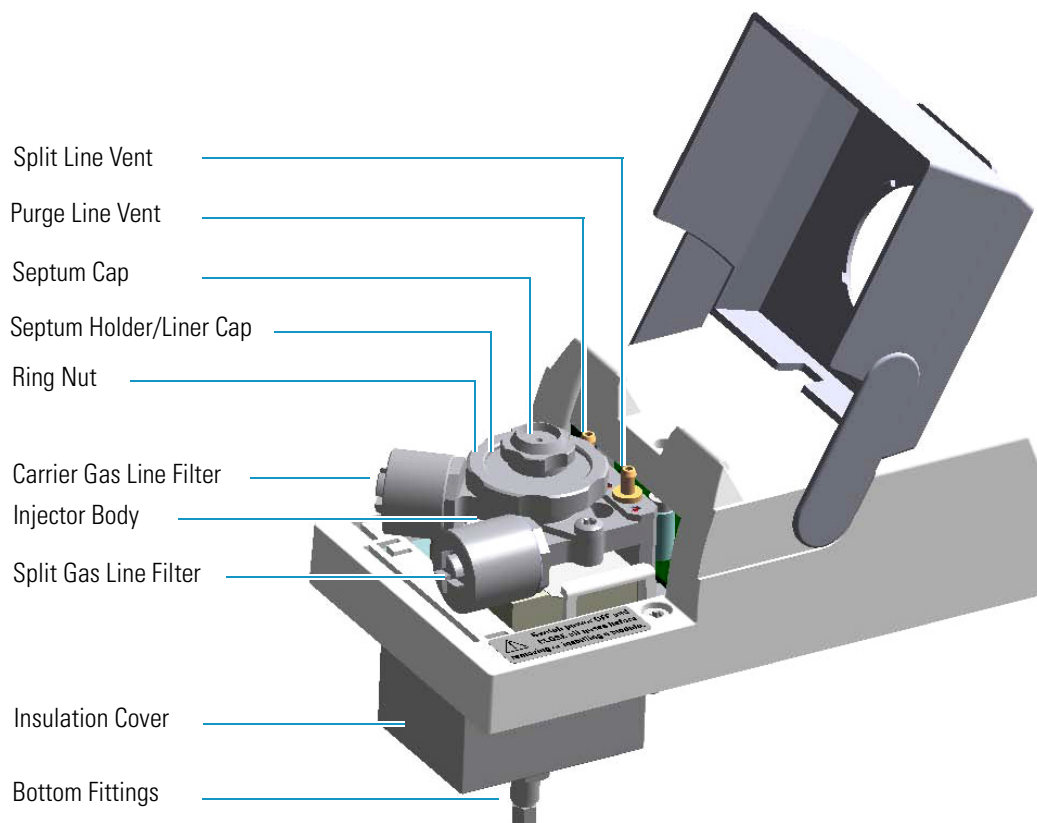

**Figure 76.** SSL Injector: Components

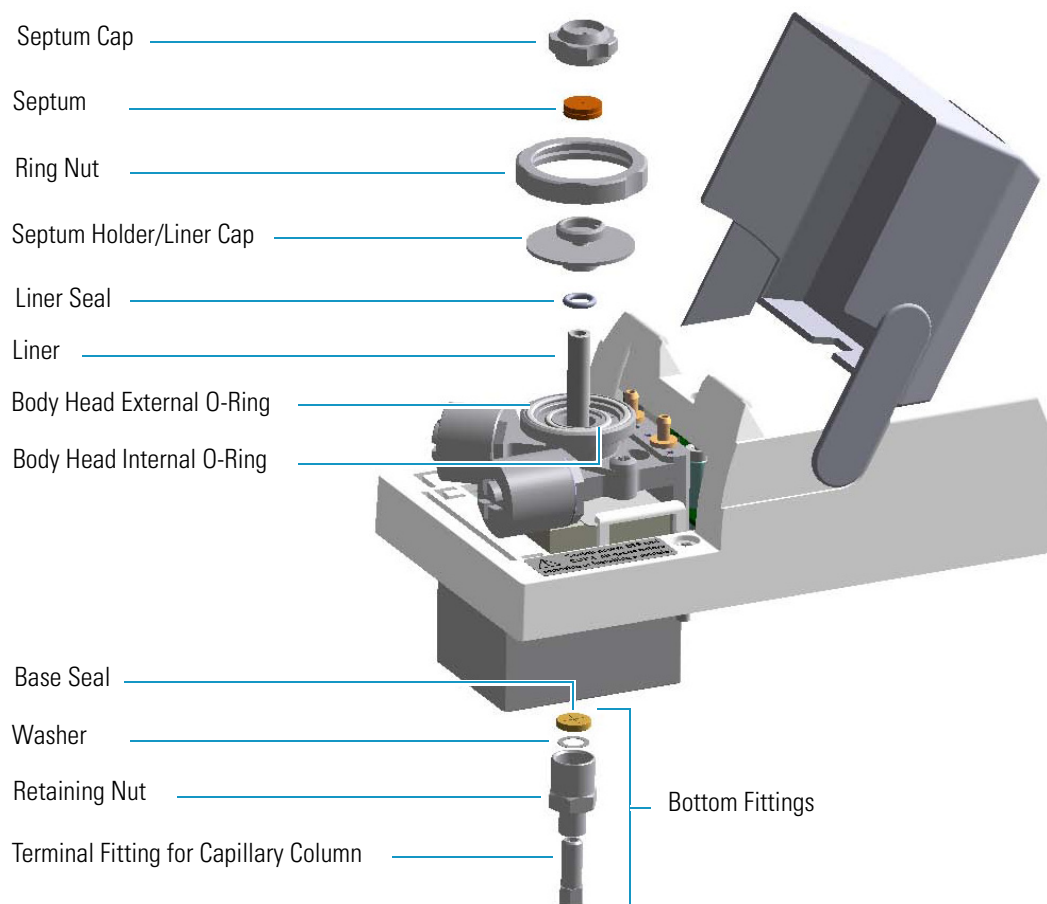

The SSL injector periodic maintenance includes:

- **Replacing the SSL septum**

Replace the septum at least after every 200 injections, or every time a problem related to septum damage, or wear occurs.

See [“Replacing the SSL Septum”](#) on [page 83](#).

- **Cleaning or replacing the SSL liner**

Check the liner for contaminants, debris, breakage, and proper installation.

The liner must be replaced depending on the number of injections performed, and the characteristics of the samples injected. Typical symptoms will indicate that the liner must be replaced. The most common is the appearance of tailing peaks in the chromatogram, particularly for polar compounds.

See [“Cleaning or Replacing the SSL Glass Liner”](#) on [page 84](#).

**Tip** It is good practice to replace the septum every time you replace the glass liner.

### 3 Performing Injectors Routine Maintenance

#### Maintaining a Split/Splitless Injector (SSL)

When replacing or removing a glass liner, it might break inside the injector. In this case the broken parts of the liner must be removed from the injector, including the glass splinters that might fall into the lower part of the vaporization chamber.

See [“Replacing a SSL Broken Liner”](#) on [page 86](#).

- **Replacing the active carbon filters on the carrier gas line and split line**

The active carbon filters must be replaced depending on the volume of solvent injected in the time.

See [“Replacing the SSL Carrier and Split Lines Filters”](#) on [page 88](#).

- **Replacing the body head o-rings**

The internal (carrier line) and external (purge line) o-rings of the body head must be replaced when in presence of leaks.

See [“Replacing the SSL Body Head O-Rings”](#) on [page 89](#).

Before maintaining the injector, read the following warning:



**WARNING** The injector fittings could be hot. Carry out all the operation at low temperature to avoid burns. Therefore, before beginning the sequence, the injector must be cooled to room temperature.



**CAUTION** When handling organic solvents, you must take precautions to avoid health hazards.

#### Materials needed to maintain a SSL injector

---

Septum

---

Tweezers

---

Glass liner

---

Liner seal (O-ring)

---

Ultrasonic cleaner

---

Mixture 1:1 methanol/acetone

---

Base Seal (if necessary)

---

Carrier gas line and/or split gas line active carbon filters

---

Body head internal o-ring

---

Body head external o-ring

---

T20 screwdriver

---




**Note** For maintaining or replacing any other component not listed in this section, see Chapter 6, “Injectors Advanced Maintenance.”

- “Baking-out Contaminants from SSL, SSLBKF, HeS-S/SL, PTV, and PTVBKF Injectors” on page 230
- “Removing/Replacing an Injector Module” on page 231
- “Cleaning the SSL Injector Body” on page 234

## Replacing the SSL Septum

### ❖ To replace the septum

**Figure 77.** SSL Injector: Septum Replacement

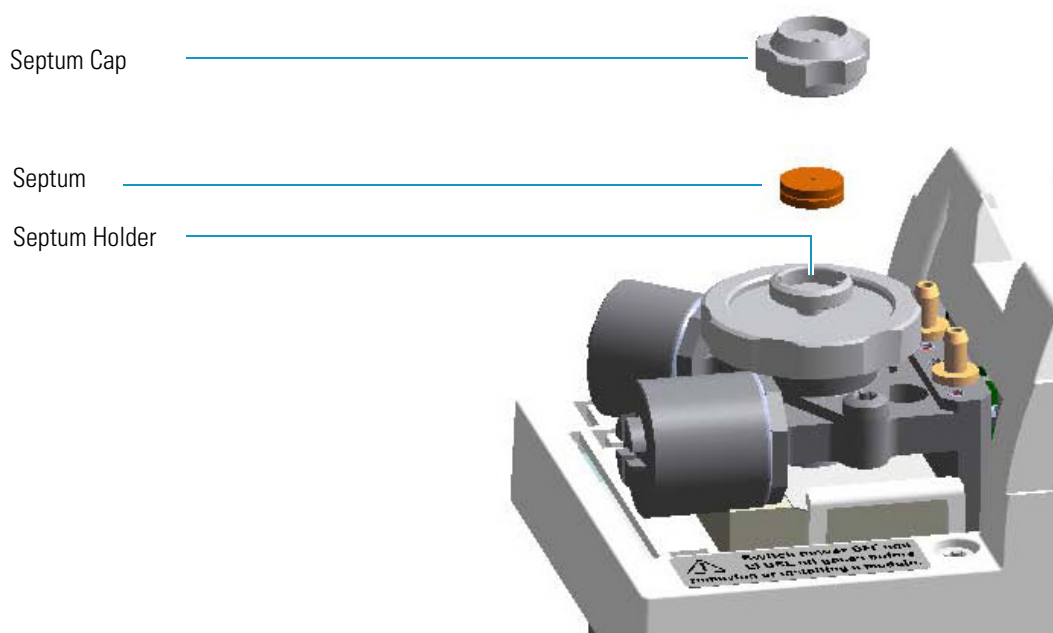

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Replace the septum.
  - a. Unscrew and remove the septum cap.

### 3 Performing Injectors Routine Maintenance

#### Maintaining a Split/Splitless Injector (SSL)

- b. Using tweezers, remove the septum from the septum holder.
- c. Avoid touching the septum with your fingers. Insert a new septum into the septum holder using tweezers.
- d. Screw and tighten the septum cap to finger-tight.




**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

7. Close the module flap cover.
8. If present, move the autosampler towards the module to restore the original alignment.
9. Turn the carrier gas on.
10. Set the normal injector, detector, and GC working conditions.

## Cleaning or Replacing the SSL Glass Liner

### ❖ To clean or replace the glass liner

**Figure 78.** SSL Injector: Glass Liner Replacement

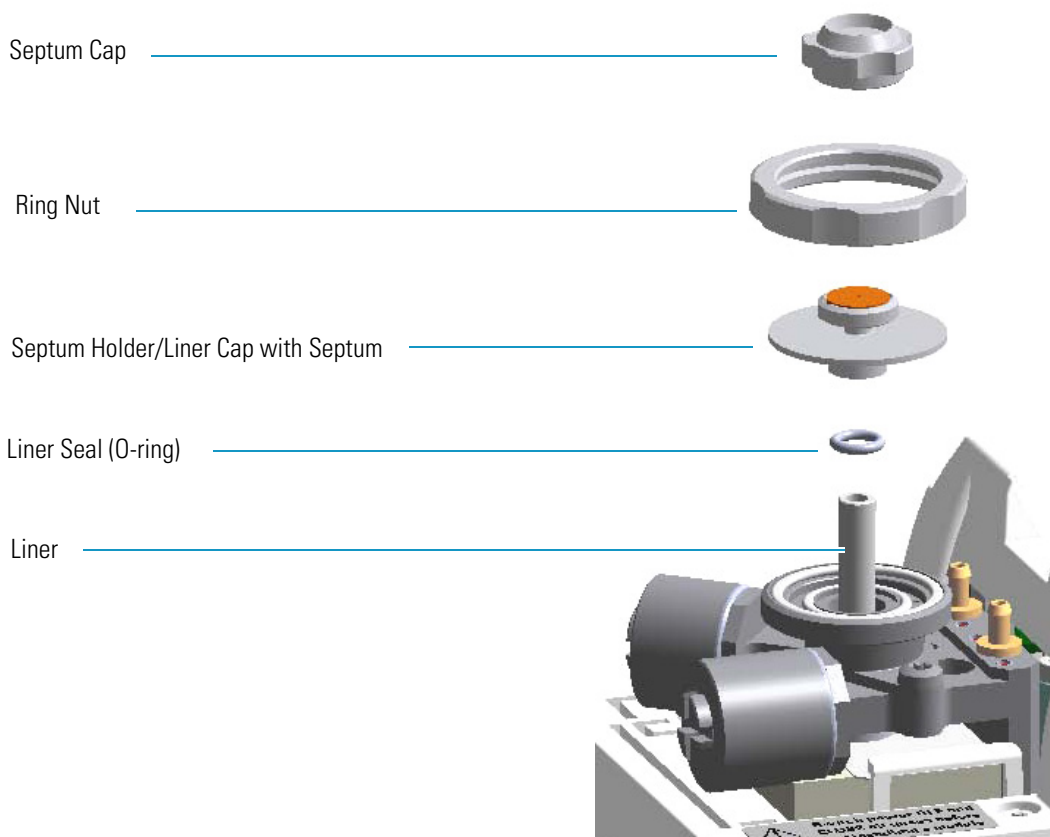

1. Put the GC in standby condition.

2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the top parts of the injector.
  - a. Unscrew the septum cap of the injector.
  - b. Unscrew the ring nut.
  - c. Remove the septum holder/liner cap with septum from the injector body head.
7. Remove the liner.
  - a. Using tweezers, remove the liner with the liner seal (o-ring) from the injector.



**CAUTION** Be careful not to break the glass liner when removing it. Glass splinters might fall into the lower part of the vaporization chamber. If the glass liner breaks, see [“Replacing a SSL Broken Liner”](#) on [page 86](#).

8. Replace or clean the liner.
  - If you are going to use a new liner, go directly to [step 10](#).
  - If you are going to clean the dirty liner, go to [step 9](#).
9. Clean the liner.
  - a. Put the liner into an ultrasonic cleaner filled with a methanol/acetone mixture (1:1).
  - b. Sonicate the liner for about half an hour.
  - c. Using tweezers, remove the liner from the bath and dry it with compressed clean air.
10. Install the liner.
  - a. Holding the new (or cleaned) liner with tweezers place a new liner seal over the liner.
  - b. Insert the liner into the injector, and push it gently towards the bottom of the injector.
11. Reinstall the top parts of the injector.
  - a. Place the septum holder/ liner cap with the septum on the body head of the injector, and fix it by screwing the ring nut.
  - b. Screw and tighten the septum cap to finger-tight.

**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

### 3 Performing Injectors Routine Maintenance


Maintaining a Split/Splitless Injector (SSL)

12. Close the module flap cover.
13. If present, move the autosampler towards the module to restore the original alignment.
14. Turn the carrier gas on.
15. Set the normal injector, detector, and GC working conditions.

## Replacing a SSL Broken Liner

### ❖ To replace a broken liner

**Figure 79.** SSL Injector: Broken Liner Replacement

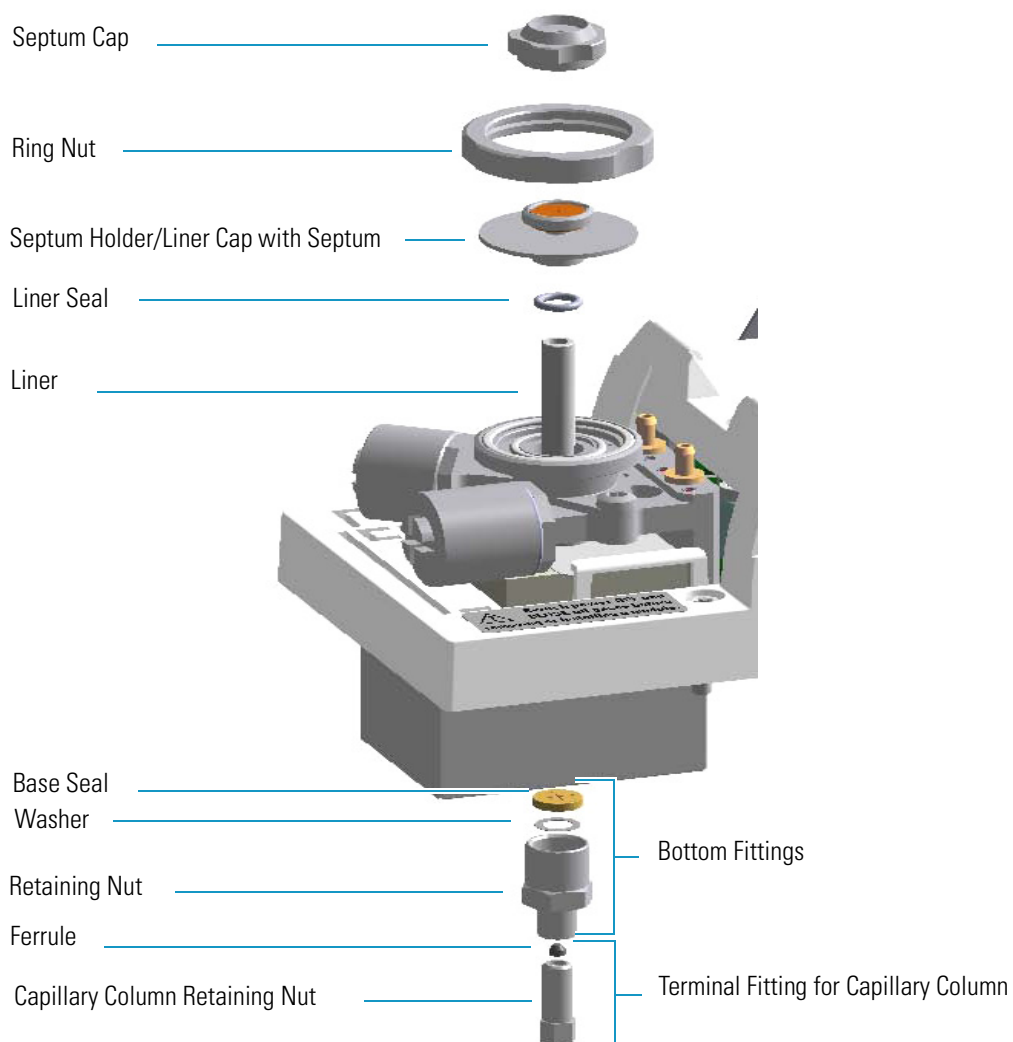

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.

3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the top parts of the injector.
  - a. Unscrew the septum cap of the injector
  - b. Unscrew the ring nut.
  - c. Remove the septum holder/liner cap with septum from the injector body head.
7. Remove the broken liner.
  - a. Using tweezers, remove the upper part of broken liner with the liner seal from the injector.
8. Remove the bottom parts of the injector.
  - a. Inside the oven, unscrew the capillary column retaining nut, then remove the analytical column with its ferrule from the bottom of the injector.
  - b. From the bottom of the injector, unscrew the retaining nut with the washer and the base seal. Glass splinters from the broken liner will fall from the injector.
  - c. With the aid of a pipe cleaner, remove the possible glass fragments from the vaporization chamber.
9. Reinstall the bottom parts of the injector.
  - a. Reinstall the analytical column.
  - b. Reinstall the retaining nut with the new washer and the base seal. If necessary replace the base seal with a new one.
10. Install a new liner.
  - a. Holding the new liner with tweezers, place a new liner seal over the liner.
  - b. Insert the liner into the injector, and push it gently towards the bottom of the injector.
11. Reinstall the top parts of the injector.
  - a. Place the septum holder/liner cap with the septum on the body head of the injector and fix them by screwing the ring nut.
  - b. Screw and tighten the injector cap to finger-tight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

### 3 Performing Injectors Routine Maintenance

Maintaining a Split/Splitless Injector (SSL)

12. Close the module flap cover.
13. If present, move the autosampler towards the module to restore the original alignment.
14. Turn the carrier gas on.
15. Set the normal injector, detector and GC working conditions.


## Replacing the SSL Carrier and Split Lines Filters



**IMPORTANT** The dimensions of the filters are different. The filter on the split gas line is bigger than the filter on the carrier gas line. Do not invert their position when you replace them. It is not necessary to replace the filters together.

### ❖ To replace the active carbon filters on carrier gas line and split line

**Figure 80.** SSL Injector: Active Filters Replacement

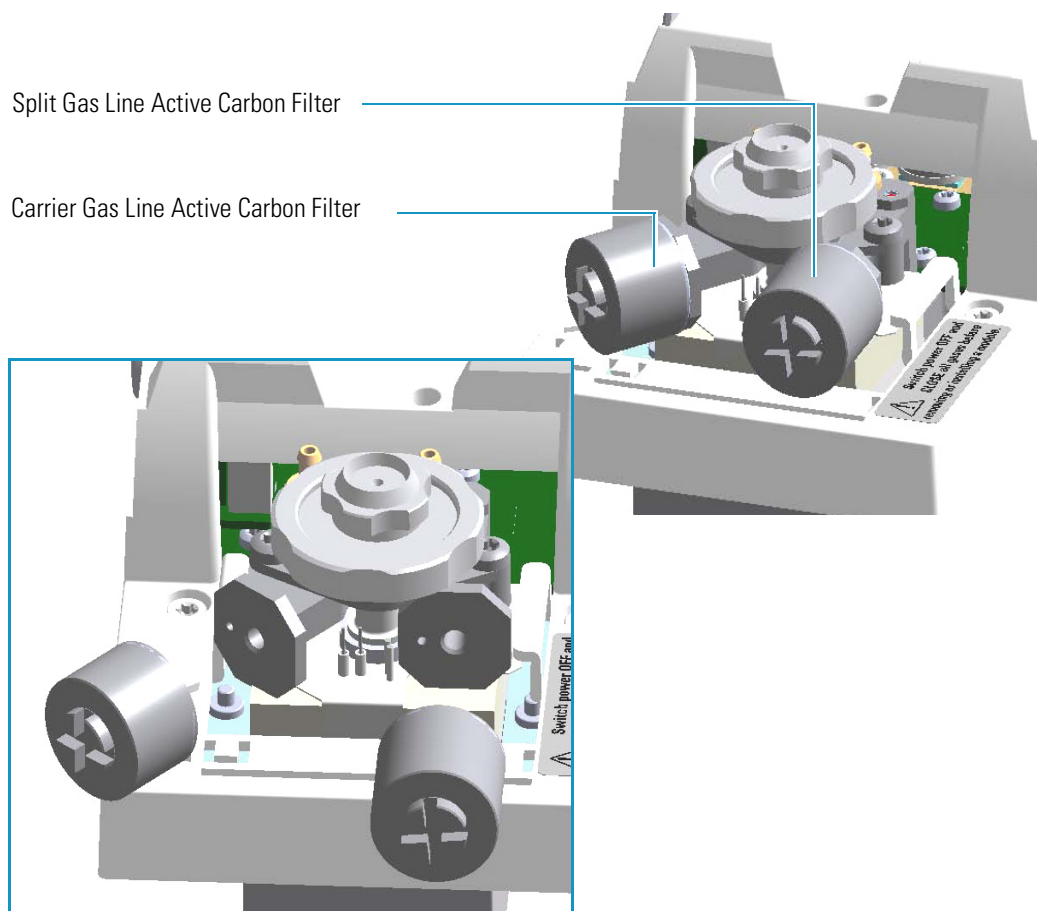

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.

3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Replace the filter.
  - a. Remove the filter to replace from its seat by turning it counter-clockwise.
  - b. Install the new filter, with o-rings, in its seat by turning it clockwise.
7. Close the module flap cover.
8. If present, move the autosampler towards the module to restore the original alignment.
9. Turn the carrier gas on.
10. Set the normal injector, detector, and GC working conditions.

## Replacing the SSL Body Head O-Rings

### ❖ To replace the body head o-rings

**Figure 81.** SSL Injector: Body Head O-Rings Replacement

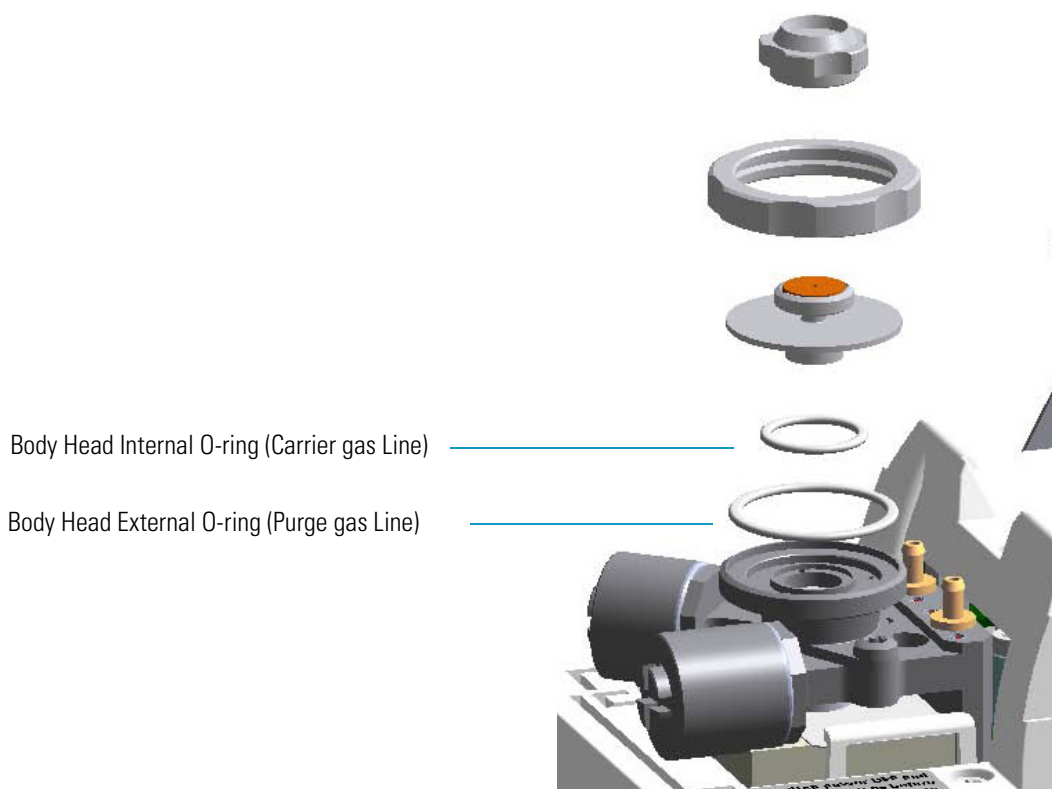

### 3 Performing Injectors Routine Maintenance

Maintaining a Split/Splitless Injector (SSL)

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the top parts of the injector.
  - a. Unscrew the septum cap of the injector.
  - b. Unscrew the ring nut.
  - c. Remove the septum holder/liner cap with septum from the injector body head.
7. Replace the head body o-rings.
  - a. Use tweezer to remove the body head internal and external o-rings, and replace them with new o-rings.
8. Reinstall the top parts of the injector.
  - a. Place the septum holder/liner cap with the septum on the body head of the injector and fix it by screwing the ring nut.
  - b. Screw and tighten the injector cap to finger-tight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

9. Close the module flap cover.
10. If present, move the autosampler towards the module to restore the original alignment.
11. Turn the carrier gas on.
12. Set the normal injector, detector, and GC working conditions.




## Maintaining a Split/Splitless Injector with Backflush (SSLBKF)

This section provides instructions for maintaining a Split/Splitless injector for Backflush (SSLBKF) applications.

The module and injector components are shown in [Figure 82](#) and [Figure 83](#).


**Figure 82.** SSLBKF Module Components

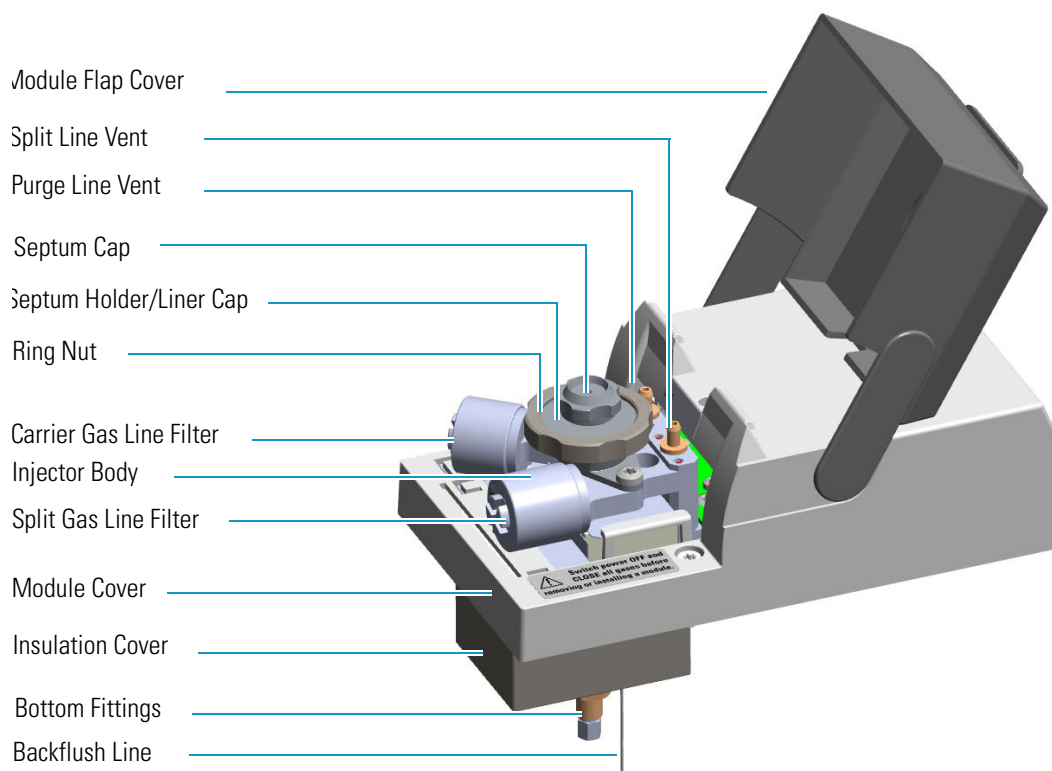

### 3 Performing Injectors Routine Maintenance

Maintaining a Split/Splitless Injector with Backflush (SSLBKF)

**Figure 83.** SSLBKF Injector Components

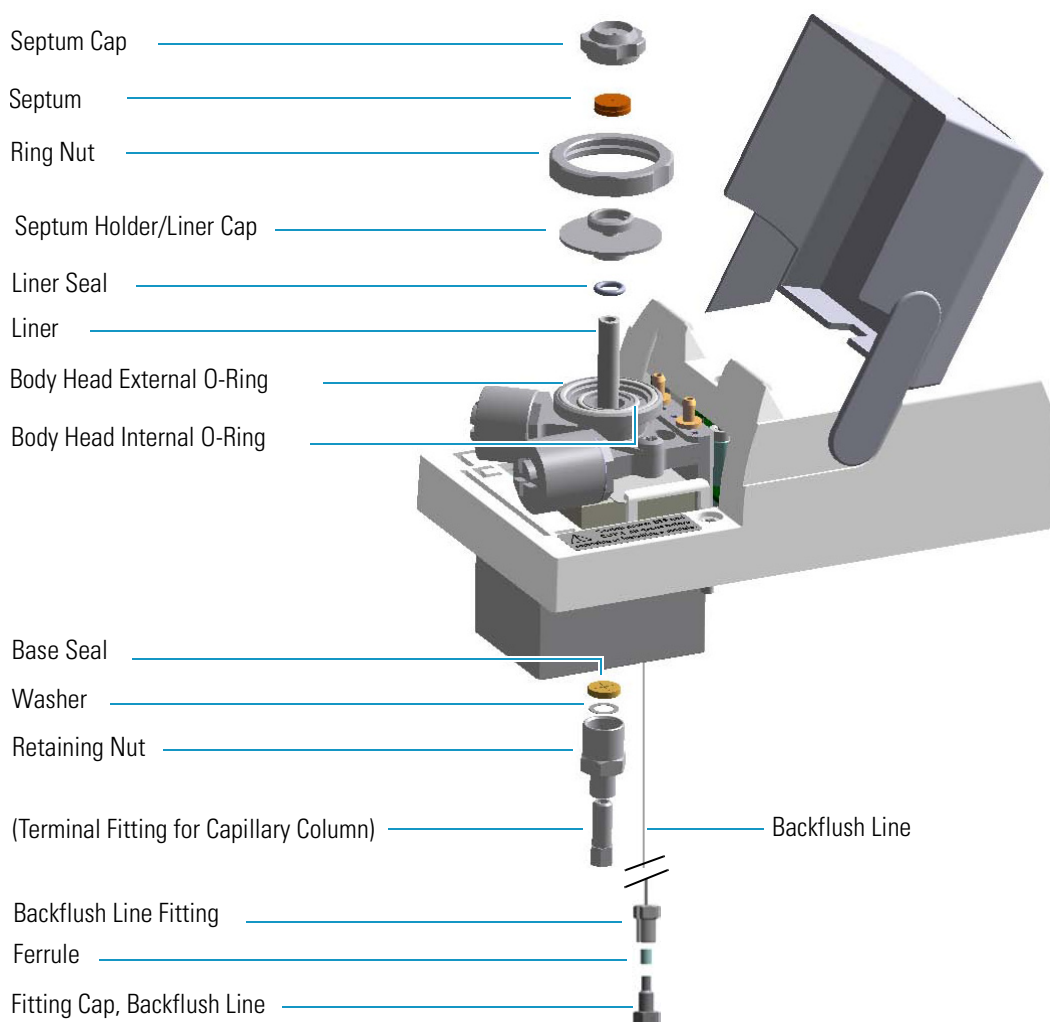

The SSLBKF injector periodic maintenance includes:

- **Replacing the SSLBKF septum**

Replace the septum at least after every 200 injections, or every time a problem related to septum damage, or wear occurs.

See [“Replacing the SSLBKF Septum”](#) on page 94.

- **Cleaning or replacing the SSLBKF liner**

Check the liner for contaminants, debris, breakage, and proper installation. The liner must be replaced depending on the number of injections performed, and the characteristics of the samples injected. Typical symptoms will indicate that the liner must be replaced. The most common is the appearance of tailing peaks in the chromatogram, particularly for polar compounds.

See [“Cleaning or Replacing the SSLBKF Glass Liner”](#) on page 95.

**Tip** It is good practice to replace the septum every time you replace the glass liner.

When replacing or removing a glass liner, it might break inside the injector. In this case the broken parts of the liner must be removed from the injector, including the glass splinters that might fall into the lower part of the vaporization chamber.

See [“Replacing a SSLBKF Broken Liner”](#) on page 97.

- **Replacing the active carbon filters on the carrier gas line and split line**

The active carbon filters must be replaced depending on the volume of solvent injected in the time.

See [“Replacing the SSLBKF Carrier and Split Lines Filters”](#) on page 99.

- **Replacing the body head o-rings**

The internal (carrier line) and external (purge line) o-rings of the body head must be replaced when in presence of leaks.

See [“Replacing the SSLBKF Body Head O-Rings”](#) on page 100.



**WARNING** The injector fittings could be hot. Carry out all the operations at low temperature to avoid burns. Therefore, before beginning the sequence, the injector must be cooled to room temperature.



**CAUTION** When handling organic solvents you must take precautions to avoid health hazards.

#### **Materials needed to maintain a SSLBKF injector**

Septum

Tweezers

Glass liner

Liner seal (O-ring)

Ultrasonic cleaner

Mixture 1:1 methanol/acetone

Base Seal (if necessary)

Carrier gas line and/or split gas line active carbon filters

Body head internal o-ring

Body head external o-ring

T20 screwdriver

### 3 Performing Injectors Routine Maintenance

Maintaining a Split/Splitless Injector with Backflush (SSLBKF)


**Note** For maintaining or replacing any other component not listed in this section, see Chapter 6, “Injectors Advanced Maintenance.”

- “Baking-out Contaminants from SSL, SSLBKF, HeS-S/SL, PTV, and PTVBKF Injectors” on page 230
- “Removing/Replacing an Injector Module” on page 231
- “Cleaning the SSLBKF Injector Body” on page 238

## Replacing the SSLBKF Septum

### ❖ To replace the septum

**Figure 84.** SSLBKF Injector: Septum Replacement

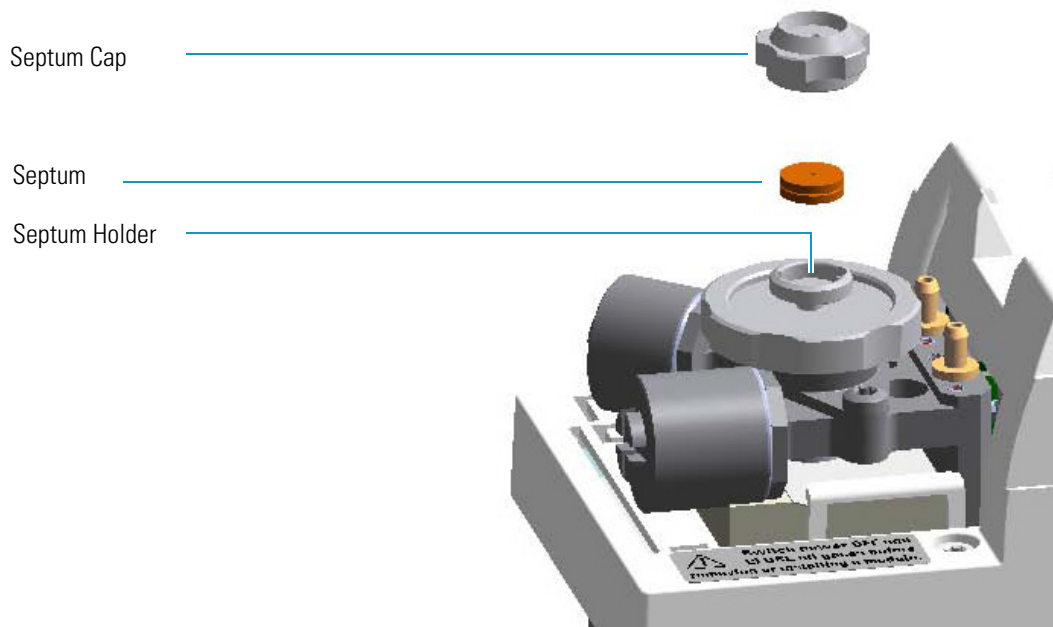

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Replace the septum.
  - a. Unscrew and remove the septum cap.

- b. Using tweezers, remove the septum from the septum holder.
  - c. Avoid touching the septum with your fingers. Insert a new septum into the septum holder using tweezers.
  - d. Screw and tighten the septum cap to finger-tight.
7. Close the module flap cover.




**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

8. If present, move the autosampler towards the module to restore the original alignment.
9. Turn the carrier gas on.
10. Set the normal injector, detector, and GC working conditions.

## Cleaning or Replacing the SSLBKF Glass Liner

### ❖ To clean or replace the glass liner

**Figure 85.** SSLBKF Injector: Glass Liner Replacement

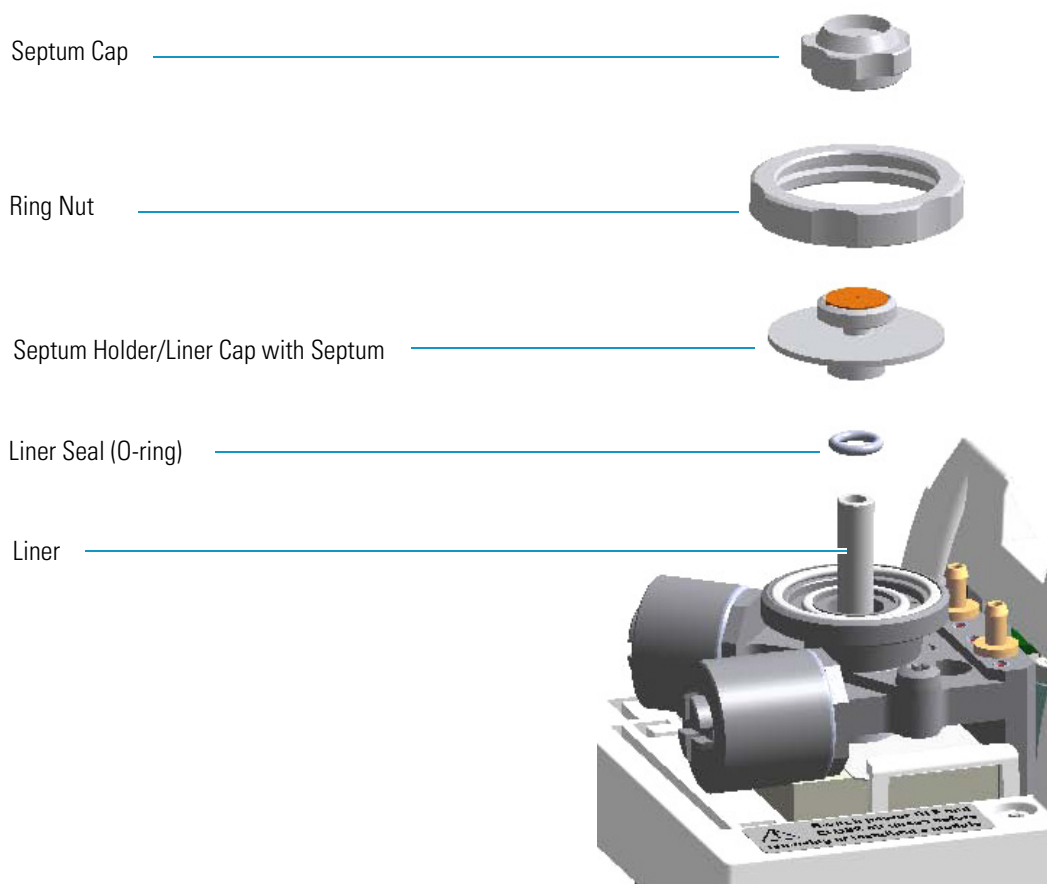

### 3 Performing Injectors Routine Maintenance

Maintaining a Split/Splitless Injector with Backflush (SSLBKF)

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the top parts of the injector.
  - a. Unscrew the septum cap of the injector.
  - b. Unscrew the ring nut.
  - c. Remove the septum holder/liner cap with septum from the injector body head.
7. Remove the liner.
  - a. Use tweezers to remove the liner with the liner seal (o-ring) from the injector.



**CAUTION** Be careful not to break the glass liner when removing it. Glass splinters might fall into the lower part of the vaporization chamber. If the glass liner breaks, see [“Replacing a SSL Broken Liner”](#) on page 86.

8. Replace or clean the liner.
  - If you are going to use a new liner, go directly to [step 10](#).
  - If you are going to clean the dirty liner, go to [step 9](#).
9. Clean the liner.
  - a. Put the liner into an ultrasonic cleaner filled with a methanol/acetone mixture (1:1).
  - b. Sonicate the liner for about half an hour.
  - c. Using tweezers, remove the liner from the bath, and dry it with compressed clean air.
10. Install the liner.
  - a. Holding the new (or cleaned) liner with tweezers place a new liner seal over the liner.
  - b. Insert the liner into the injector and push it gently towards the bottom of the injector.
11. Reinstall the top parts of the injector.
  - a. Place the septum holder/ liner cap with the septum on the body head of the injector and fix it by screwing the ring nut.
  - b. Screw and tighten the septum cap to finger-tight.




**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

12. Close the module flap cover.
13. If present, move the autosampler towards the module to restore the original alignment.
14. Turn the carrier gas on.
15. Set the normal injector, detector, and GC working conditions.

## Replacing a SSLBKF Broken Liner

### ❖ To replace a broken liner

**Figure 86.** SSLBKF Injector: Broken Liner Replacement

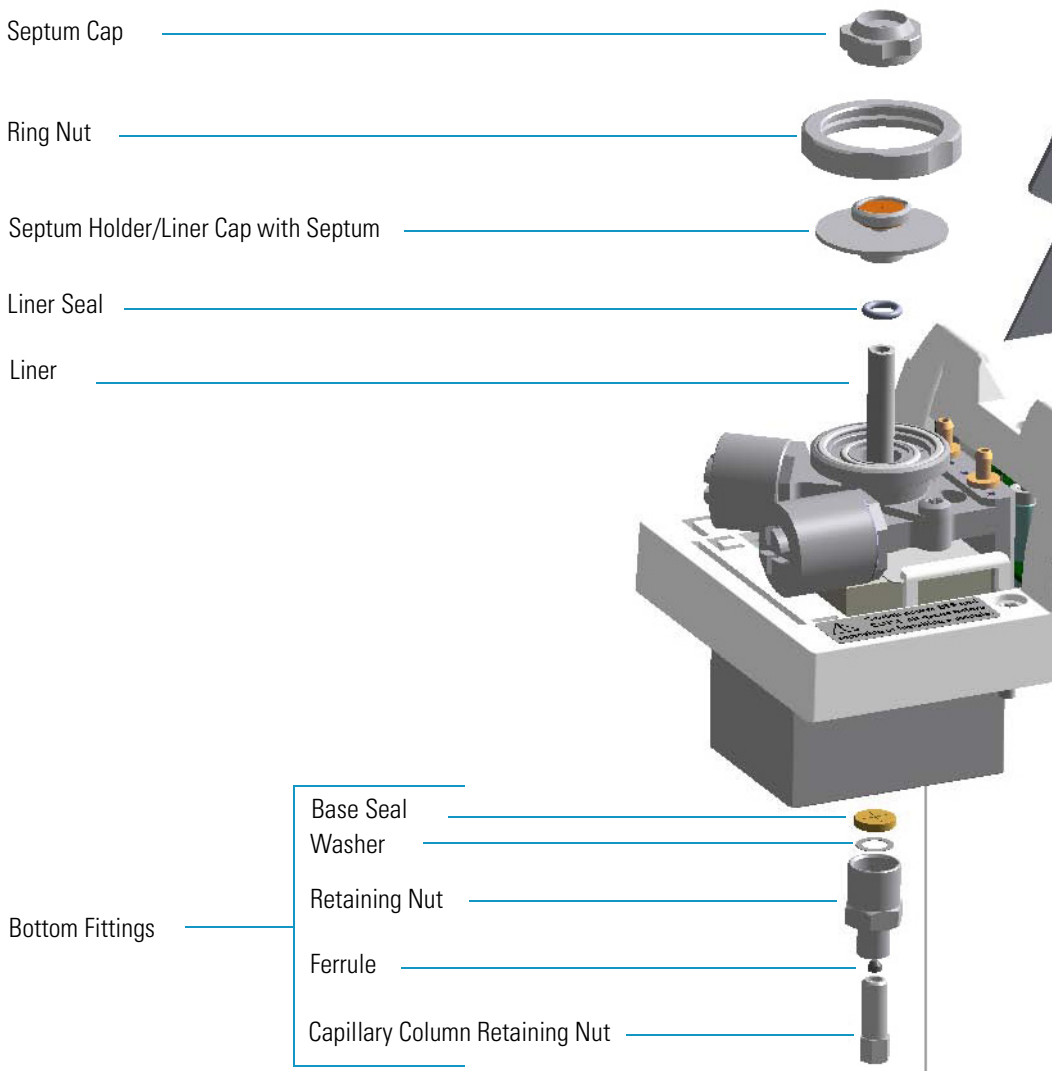

### 3 Performing Injectors Routine Maintenance

Maintaining a Split/Splitless Injector with Backflush (SSLBKF)

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the top parts of the injector.
  - a. Unscrew the septum cap of the injector.
  - b. Unscrew the ring nut.
  - c. Remove the septum holder/liner cap with septum from the injector body head.
7. Remove the broken liner.
  - a. Use tweezers to remove the upper part of broken liner with the liner seal from the injector.
8. Remove the bottom parts of the injector.
  - a. Inside the oven, unscrew the capillary column retaining nut, then remove the analytical column with its ferrule from the bottom of the injector.
  - b. From the bottom of the injector, unscrew the retaining nut with the washer and the base seal. Glass splinters from the broken liner will fall from the injector.
  - c. With the aid of a pipe cleaner, remove the possible glass fragments from the vaporization chamber.
9. Reinstall the bottom parts of the injector.
  - a. Reinstall the analytical column.
  - b. Reinstall the retaining nut with the new washer and the base seal (if necessary replace the base seal with a new one).
10. Install a new liner.
  - a. Holding the new liner with tweezers, place a new liner seal over the liner.
  - b. Insert the liner into the injector and push it gently towards the bottom of the injector.
11. Reinstall the top parts of the injector.
  - a. Place the septum holder/liner cap with the septum on the body head of the injector and fix them by screwing the ring nut.
  - b. Screw and tighten the injector cap to finger-tight





**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

12. Close the module flap cover.
13. If present, move the autosampler towards the module to restore the original alignment.
14. Turn the carrier gas on.
15. Set the normal injector, detector, and GC working conditions.


## Replacing the SSLBKF Carrier and Split Lines Filters



**IMPORTANT** The dimensions of the filters are different. The filter on the split gas line is bigger than the filter on the carrier gas line. Do not invert their position when you replace them. It is not necessary to replace the filters together.

### ❖ To replace the active carbon filters on carrier gas line and split line

**Figure 87.** SSLBKF Injector: Active Filters Replacement

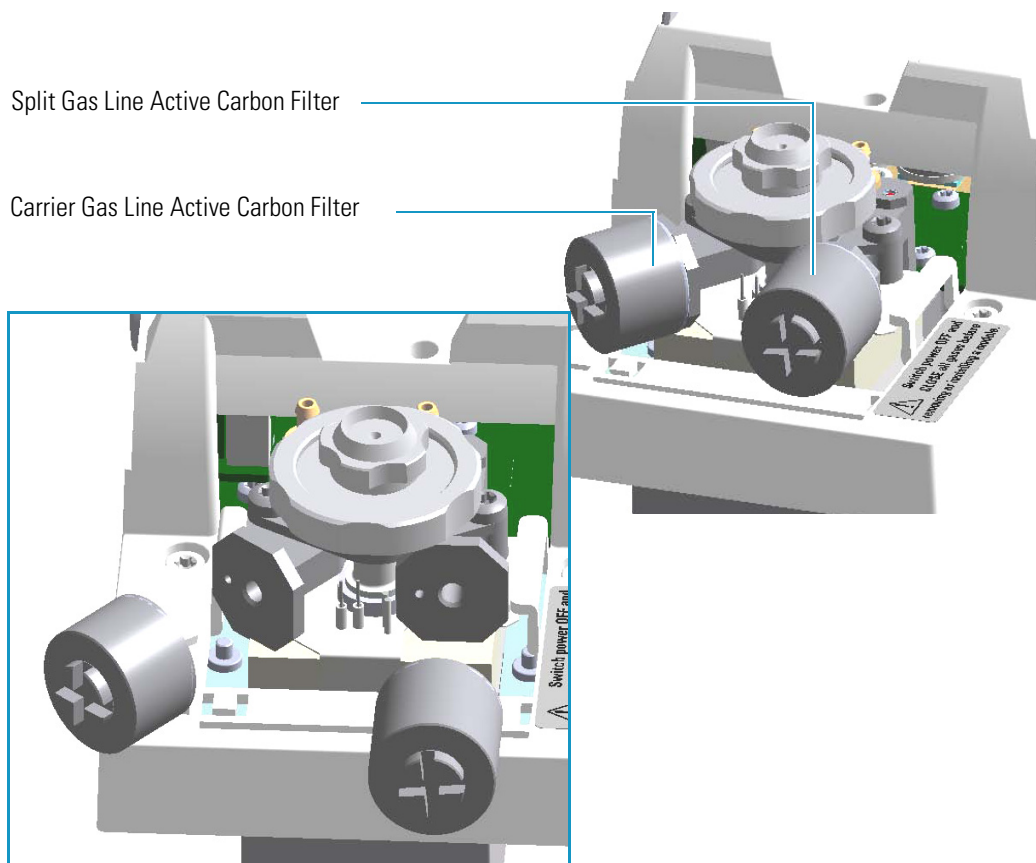

1. Put the GC in standby condition.

### 3 Performing Injectors Routine Maintenance

Maintaining a Split/Splitless Injector with Backflush (SSLBKF)

2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Replace the filter.
  - a. Remove the filter to replace from its seat by turning it counter-clockwise.
  - b. Install the new filter, with o-rings, in its seat by turning it clockwise.
7. Close the module flap cover.
8. If present, move the autosampler towards the module to restore the original alignment.
9. Turn the carrier gas on.
10. Set the normal injector, detector, and GC working conditions.

## Replacing the SSLBKF Body Head O-Rings

### ❖ To replace the body head o-rings

**Figure 88.** SSLBKF Injector: Body Head O-Rings Replacement

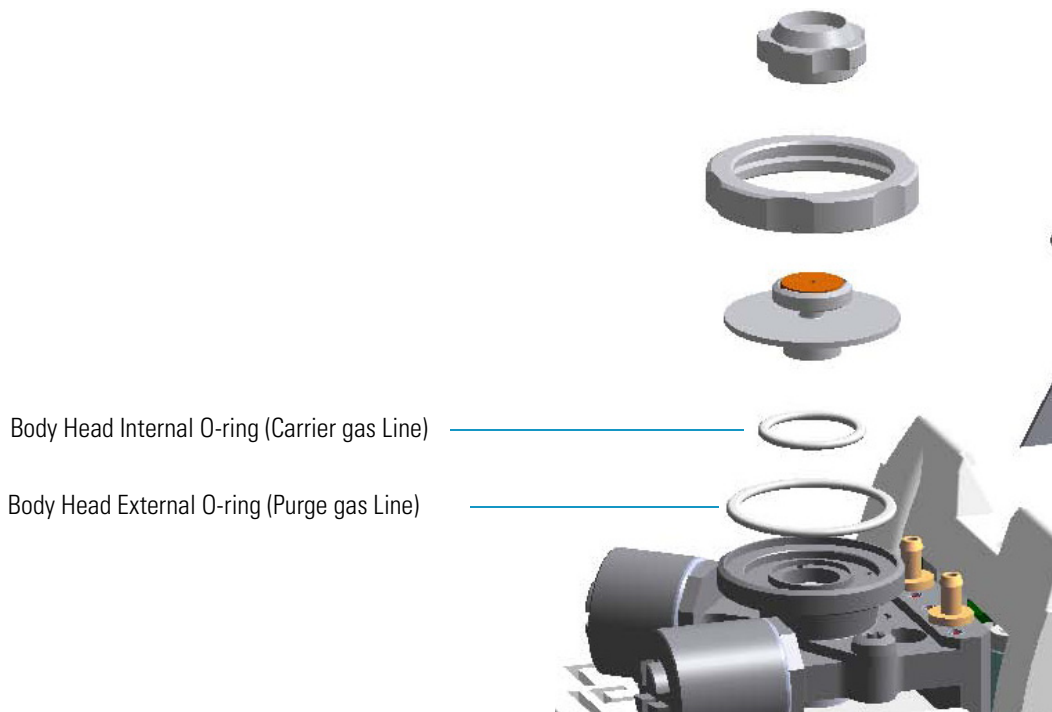

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the top parts of the injector.
  - a. Unscrew the septum cap of the injector.
  - b. Unscrew the ring nut.
  - c. Remove the septum holder/liner cap with septum from the injector body head.
7. Replace the head body o-rings.
  - a. Use tweezers to remove the body head internal and external o-rings, and replace them with new o-rings.
8. Reinstall the top parts of the injector.
  - a. Place the septum holder/liner cap with the septum on the body head of the injector and fix it by screwing the ring nut.
  - b. Screw and tighten the injector cap to finger-tight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.


9. Close the module flap cover.
10. If present, move the autosampler towards the module to restore the original alignment.
11. Turn the carrier gas on.
12. Set the normal injector, detector, and GC working conditions.

## Maintaining a Gas Sampling Valve Injector (GSV)


This section provides instructions for maintaining a Gas Sampling Valve injector (GSV) module.

The module and injector components are shown in [Figure 89](#).

**Figure 89.** GSV Module Components

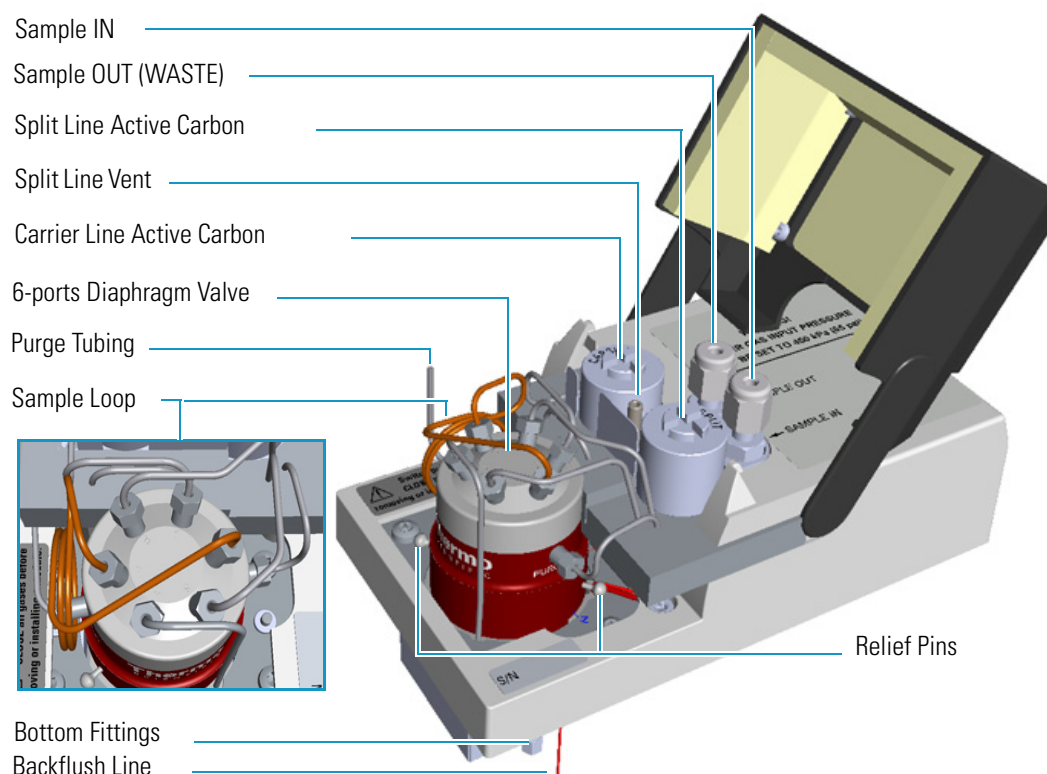

**IMPORTANT** The valve is delivered with relief pins. These pins **MUST BE** removed from the valve body before working with the valve.

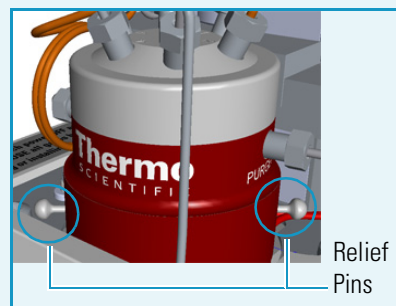

Once actuation tubing has been set up, and pressure has been adjusted to 450 kPa (65 psig), actuate the valve and remove the pins. You may need to increase the pressure a little to ease hand removal of the pins.

Remember, make sure to properly readjust the actuating operating pressure after removing the pins, if you had to change it. Keep these pins in a safe place. You may want to re-use them for valve maintenance. It is good practice to re-install the relief

pins if the valve is not used for a long time.

For details refer to the instructions provided by the valve manufacturer.

The Gas Sampling Valve periodic maintenance includes:

- **Connecting the sample In and Out lines**

See [Connecting the Sample In and Out Lines](#).

- **Replacing the active carbon filters on the carrier gas line and split line**

The active carbon filters must be replaced depending on the volume of solvent injected in the time.

See [Replacing the Carrier and Split Lines Filters](#).

- **Replacing the sample loop**

Replace the sample loop when an upper volume of sample than the volume of standard loop is required. Choose among the loops available.

See [Replacing the Sample Loop](#).

- **Replacing the gas sampling valve diaphragm**


Follow the instructions provided by the manufacturer of the valve.

## Connecting the Sample In and Out Lines

❖ **To connect the sample In and Out lines**

1. If not already done, remove the protective caps from the Sample In and Sample Out fittings. See [Figure 90](#).

**Figure 90.** Protective Caps

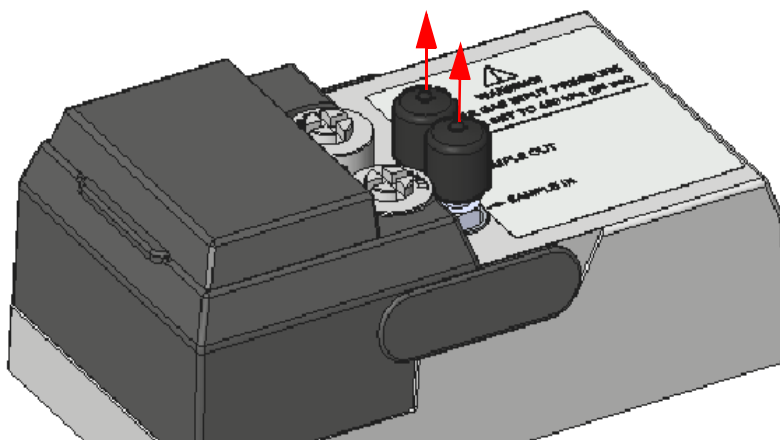


1. By using the proper 1/8-in. tubing, nut and ferrule, connect the inlet sample line to the **Sample In** port on the GSV module. See [Figure 91](#).

### 3 Performing Injectors Routine Maintenance


Maintaining a Gas Sampling Valve Injector (GSV)

**Figure 91.** GSV Sample In Line Connection

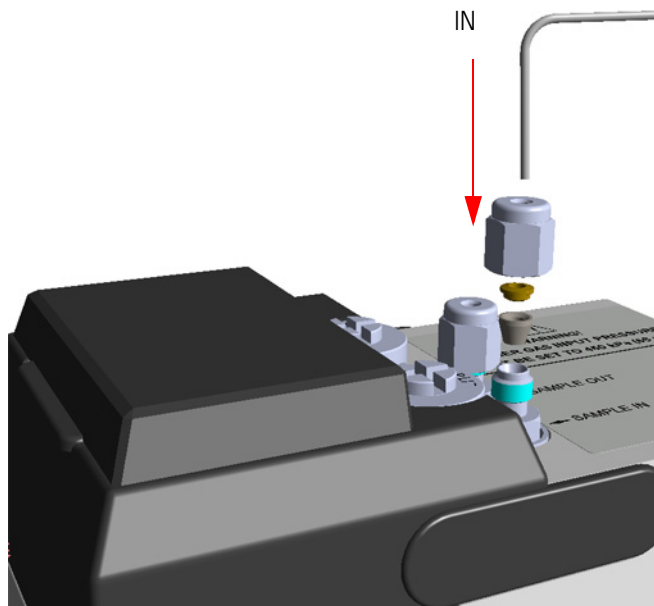


2. By using the proper 1/8-in. tubing, nut and ferrule, connect the **Sample OUT** port on the GSV module to an exhaust system. See [Figure 92](#).

**Figure 92.** GSV Sample Out Line Connection

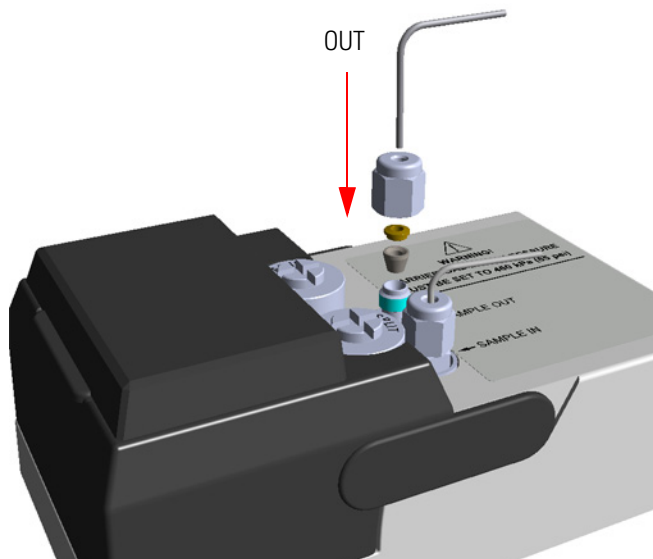

**Figure 93.** GSV Sample In and Out Lines Connection

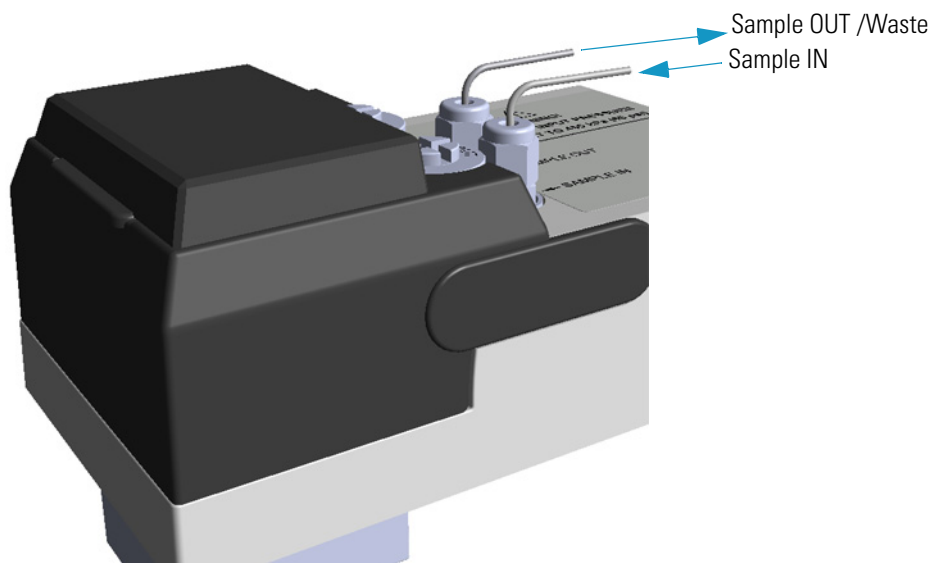

## Replacing the Carrier and Split Lines Filters



**IMPORTANT** The dimensions of the filters are different. The filter on the split gas line is bigger than the filter on the carrier gas line. Do not invert their position when you replace them. It is not necessary to replace the filters together.

### ❖ To replace the active carbon filters on carrier gas line and split line

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Replace the filter.
  - a. Remove the filter to replace from its seat by turning it counter-clockwise.
  - b. Install the new filter, with o-rings, in its seat by turning it clockwise.



### 3 Performing Injectors Routine Maintenance

Maintaining a Gas Sampling Valve Injector (GSV)

**Figure 94.** GSV Module Active Filters Replacement

Split Gas Line Active Carbon Filter

Carrier Gas Line Active Carbon Filter

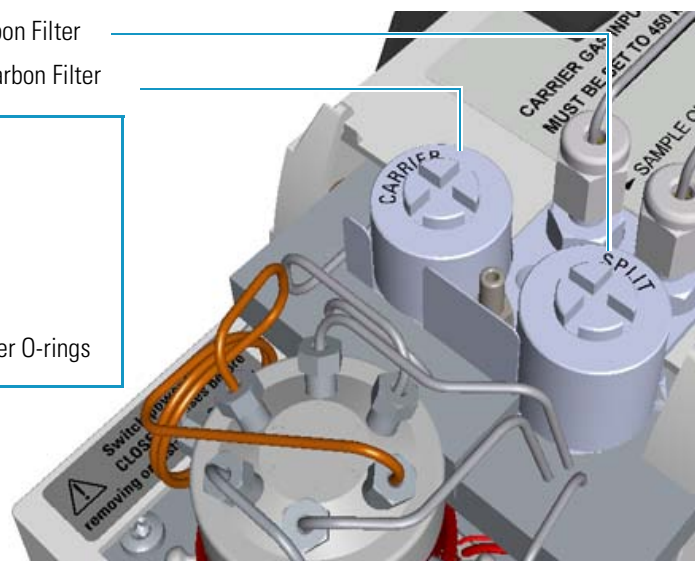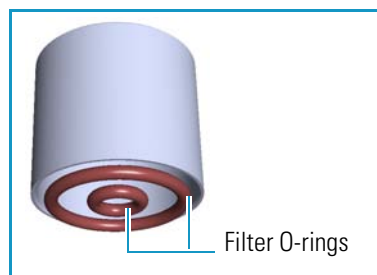

7. Close the module flap cover.
8. If present, move the autosampler towards the module to restore the original alignment.
9. Turn the carrier gas on.
10. Set the normal injector, detector, and GC working conditions.

## Replacing the Sample Loop

### ❖ To replace the sampling loop

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Replace the sampling loop.
  - a. Use a 3/16-in. wrench to loosen the two nuts connecting the sample loop to the six-port valve.
  - b. Remove the loop off the valve.
  - c. Mount the new loop over the valve and tighten the nuts to the six-port valve.



**Figure 95.** Example of Sample Loop Replacement

Sample Loop

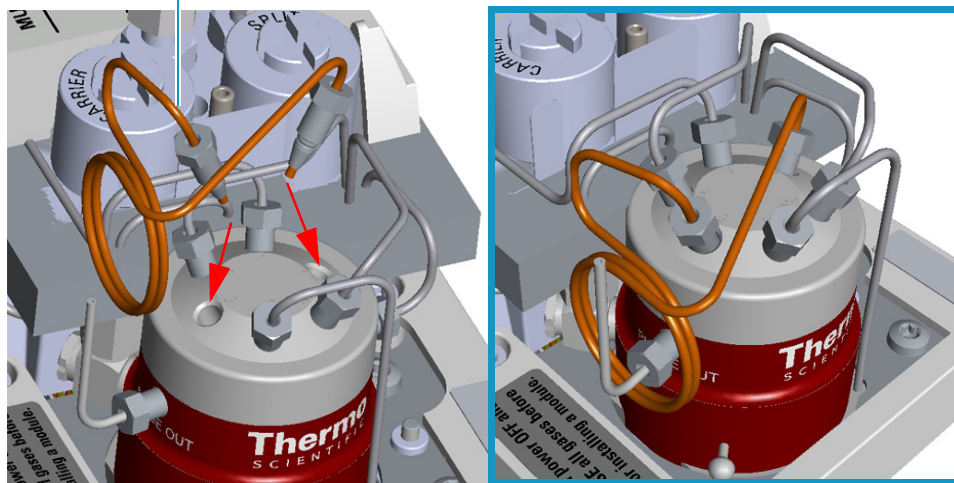

7. Close the module flap cover.
8. If present, move the autosampler towards the module to restore the original alignment.
9. Turn the carrier gas on.
10. Set the normal injector, detector, and GC working conditions.

### 3 Performing Injectors Routine Maintenance


Maintaining an Instant Connect Helium Saver Injector (HeS-S/SL)

## Maintaining an Instant Connect Helium Saver Injector (HeS-S/SL)


This section provides instructions for maintaining an Instant Connect Helium Saver Injector Module (HeS-S/SL).

The module and injector components are shown in [Figure 96](#) and [Figure 97](#).

**Figure 96.** HeS-S/SL Module Components

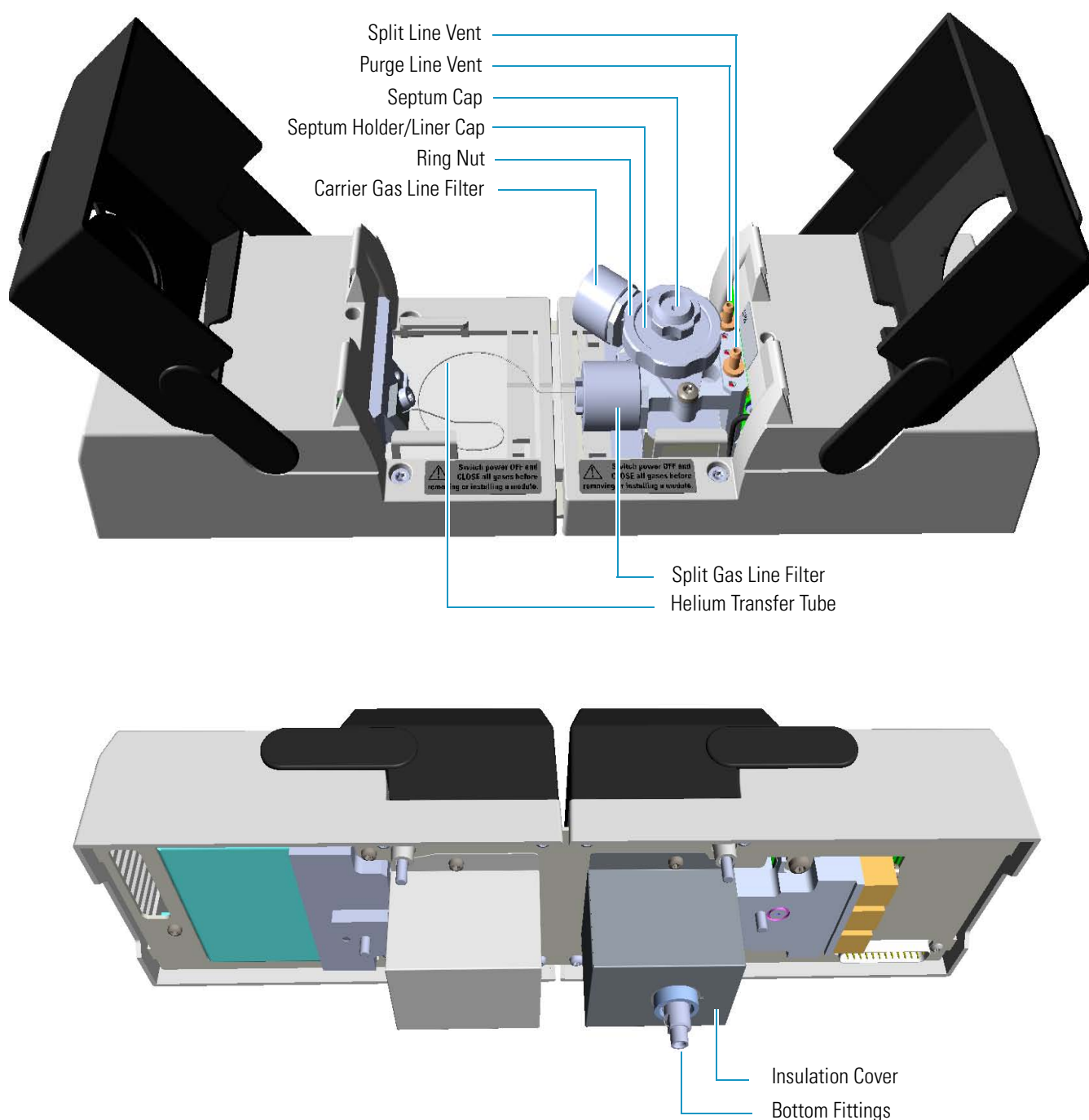

**Figure 97.** HeS-S/SL Module Components

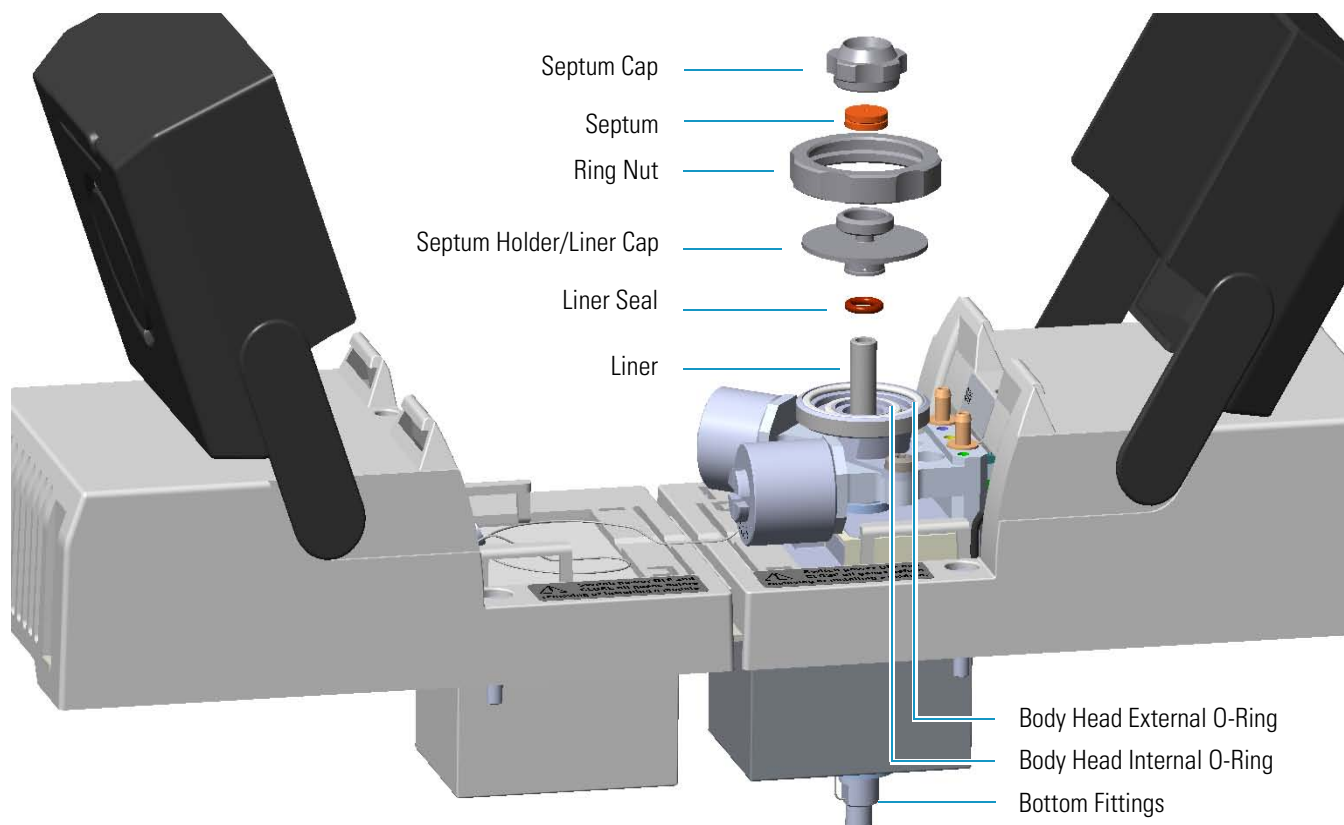

Maintaining the Instant Connect Helium Saver Injector Module is largely the same as a conventional SSL inlet.

One advantage of the Instant Connect Helium Saver Injector Module is that routine septum and liner changes can be accomplished without cooling the MS transfer line or ion source. This is particularly desirable since cooling down and re-establishing stable MS temperatures takes much longer than cooling down and re-establishing the inlet temperature.

The Instant Connect Helium Saver Injector Module periodic maintenance includes:

- **Replacing the septum**

The septum needs to be changed intermittently to prevent leakage. Replace the septum at least after every 200 injections, or every time a problem related to septum damage, or wear occurs.

See [“Replacing the Septum”](#) on [page 111](#).

- **Cleaning or replacing the liner**

Injection port liner needs to be replaced or cleaned as it becomes dirty.

See [“Cleaning or Replacing the Glass Liner”](#) on [page 112](#).

**Tip** It is good practice to replace the septum every time you replace the glass liner.

### 3 Performing Injectors Routine Maintenance

Maintaining an Instant Connect Helium Saver Injector (HeS-S/SL)

- **Replacing the active carbon filters on the carrier gas line and split line**

The active carbon filters must be replaced depending on the volume of solvent injected in the time.

See “Replacing the Carrier and Split Lines Filters” on page 114.

- **Replacing the body head o-rings**

The internal (carrier line) and external (purge line) o-rings of the body head must be replaced when there are leaks present.

See “Replacing the Body Head O-Rings” on page 115.

Before maintaining the injector, read the following warning:



**WARNING** The injector fittings could be hot. Carry out all the operation at low temperature to avoid burns. Therefore, before beginning the sequence, the injector must be cooled to room temperature.



**CAUTION** When handling organic solvents, you must take precautions to avoid health hazards.

#### Materials needed to maintain a HeS-S/SL

Septum

Tweezers

Glass liner

Liner seal (O-ring)

Ultrasonic cleaner

Mixture 1:1 methanol/acetone

Carrier gas line and/or split gas line active carbon filters

Body head internal o-ring

Body head external o-ring

T20 screwdriver


**Note** For maintaining or replacing any other component not listed in this section, see Chapter 6, “Injectors Advanced Maintenance.”

- “Baking-out Contaminants from SSL, SSLBKF, HeS-S/SL, PTV, and PTVBKF Injectors” on page 230
- “Cleaning the HeS-S/SL Injector Body” on page 242

## Replacing the Septum

### ❖ To replace the septum

**Figure 98.** HeS-S/SL Injector: Septum Replacement

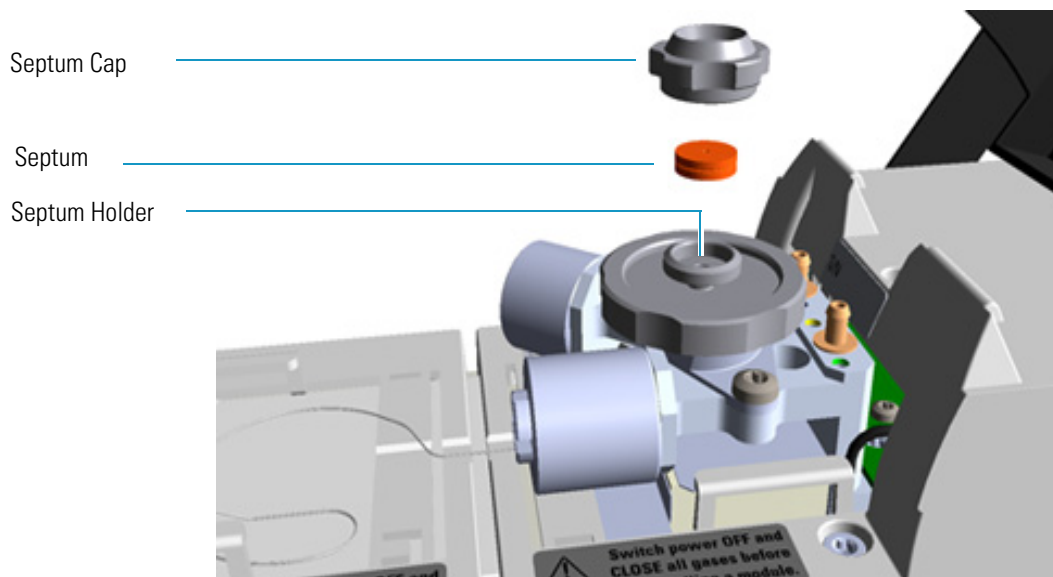

**Note** Care should be taken when performing this procedure in order to keep from damaging analytical columns.

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Ensure the **Enable conservation** field in the configuration page of the inlet reads **Yes**.
4. Check that there is at least 20 mL/min split flow exiting the split line, then turn the column flow to **Off** in the GC user interface and allow the inlet to depressurize.
5. Put the autosampler away if present.
6. Open the module flap cover.
7. Replace the septum.
  - a. Unscrew and remove the septum cap.
  - b. Using tweezers, remove the septum from the septum holder.
  - c. Avoid touching the septum with your fingers. Insert a new septum into the septum holder using tweezers.
  - d. Screw and tighten the septum cap to finger-tight.

### 3 Performing Injectors Routine Maintenance

Maintaining an Instant Connect Helium Saver Injector (HeS-S/SL)




**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

8. Close the module flap cover.
9. If present, move the autosampler towards the module to restore the original alignment.
10. Wait two or three minutes for helium to purge the inlet.
11. Turn the inlet flow back on and enable the inlet heater.
12. Set the normal injector, detector, and GC working conditions.

## Cleaning or Replacing the Glass Liner

### ❖ To clean or replace the glass liner

**Figure 99.** HeS-S/SL Injector: Glass Liner Replacement

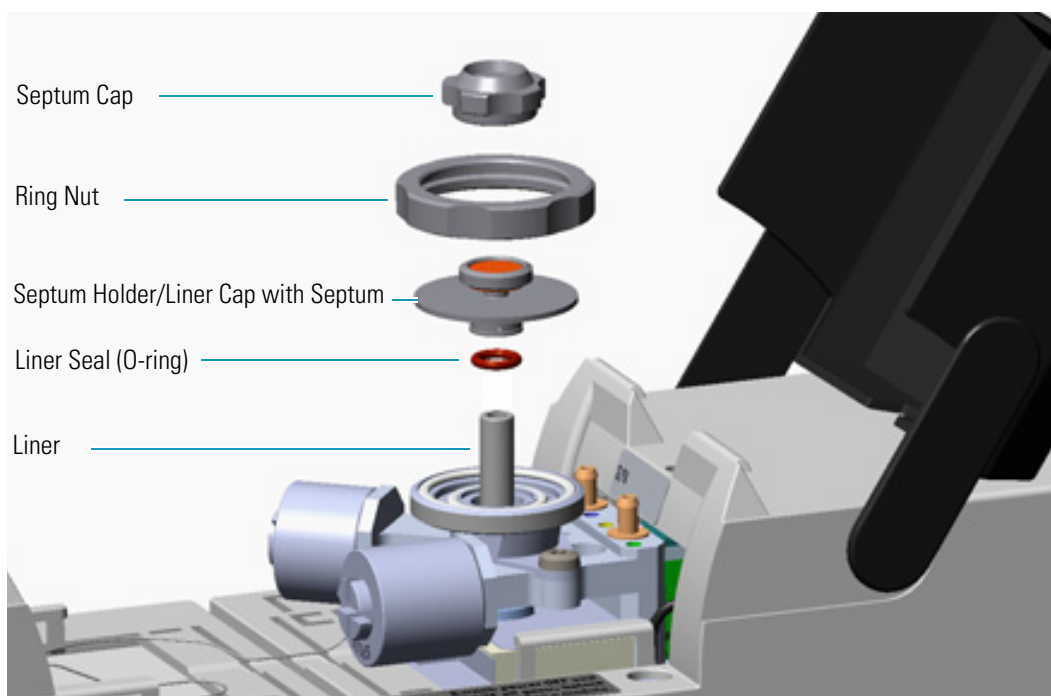

**Note** Care should be taken when performing the procedure in order to keep from damaging analytical columns.

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Ensure the **Enable conservation** field in the configuration page of the inlet reads **Yes**.
4. Check that there is at least 20 mL/min split flow exiting the split line, then turn the column flow to **Off** in the GC user interface and allow the inlet to depressurize.
5. Put the autosampler away if present.
6. Open the module flap cover.
7. Remove the top parts of the injector.
  - a. Unscrew the septum cap of the injector.
  - b. Unscrew the ring nut.
  - c. Remove the septum holder/liner cap with septum from the injector body head.
8. Remove the liner.
  - a. Using tweezers, remove the liner with the liner seal (o-ring) from the injector.
9. Replace or clean the liner.
  - If you are going to clean the dirty liner, go to [step 10](#).
  - If you are going to use a new liner, go directly to [step 11](#).
10. Clean the liner.
  - a. Put the liner into an ultrasonic cleaner filled with a methanol/acetone mixture (1:1).
  - b. Sonicate the liner for about half an hour.
  - c. Using tweezers, remove the liner from the bath and dry it with compressed clean air.
11. Install the liner.
  - a. Hold the new (or cleaned) liner with tweezers place a new liner seal over the liner.
  - b. Insert the liner into the injector, and push it gently towards the bottom of the injector.
12. Reinstall the top parts of the injector.
  - a. Place the septum holder/liner cap with the septum on the body head of the injector, and fix it by screwing the ring nut.
  - b. Screw and tighten the septum cap to finger-tight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

13. Close the module flap cover.
14. If present, move the autosampler towards the module to restore the original alignment.
15. Wait two or three minutes for helium to purge the inlet.

### 3 Performing Injectors Routine Maintenance

Maintaining an Instant Connect Helium Saver Injector (HeS-S/SL)

16. Turn the inlet flow back on and enable the inlet heater.
17. Set the normal injector, detector, and GC working conditions.

## Replacing the Carrier and Split Lines Filters




**CAUTION** The dimensions of the filters are different. The filter on the split gas line is bigger than the filter on the carrier gas line. Do not invert their position when you replace them. It is not necessary to replace the filters together.

### ❖ To replace the active carbon filters on carrier gas line and split line

**Figure 100.** HeS-S/SL Injector: Active Filters Replacement

Split Gas Line Active Carbon Filter  
Carrier Gas Line Active Carbon Filter

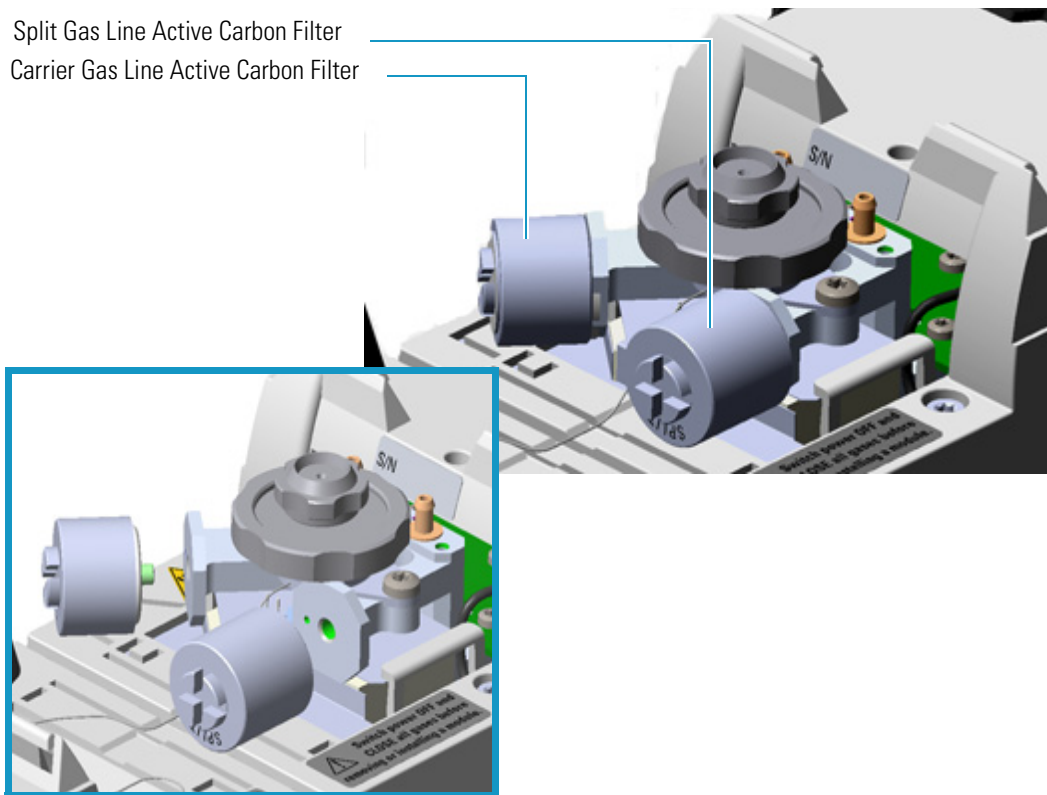

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Ensure the **Enable conservation** field in the configuration page of the inlet reads **Yes**.
4. Check that there is at least 20 mL/min split flow exiting the split line, then turn the column flow to **Off** in the GC user interface and allow the inlet to depressurize.




5. Put the autosampler away if present.
6. Open the module flap cover.
7. Replace the filter.
  - a. Remove the filter to replace from its seat by turning it counter-clockwise.
  - b. Install the new filter, with o-rings, in its seat by turning it clockwise.
8. Close the module flap cover.
9. If present, move the autosampler towards the module to restore the original alignment.
10. Wait two or three minutes for helium to purge the inlet.
11. Turn the inlet flow back on and enable the inlet heater.
12. Set the normal injector, detector, and GC working conditions.

## Replacing the Body Head O-Rings

### ❖ To replace the body head o-rings

**Figure 101.** HeS-S/SL Injector: Body Head O-Rings Replacement

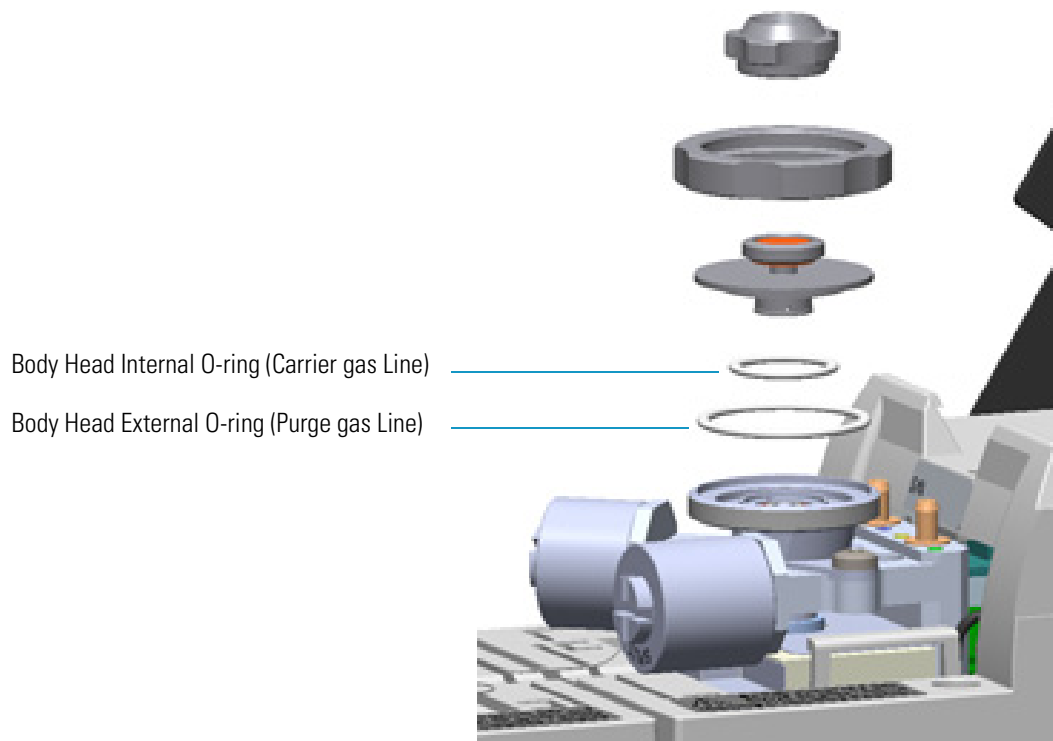

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.

### 3 Performing Injectors Routine Maintenance

Maintaining an Instant Connect Helium Saver Injector (HeS-S/SL)

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Ensure the **Enable conservation** field in the configuration page of the inlet reads **Yes**.
4. Check that there is at least 20 mL/min split flow exiting the split line, then turn the column flow to **Off** in the GC user interface and allow the inlet to depressurize.
5. Put the autosampler away if present.
6. Open the module flap cover.
7. Remove the top parts of the injector.
  - a. Unscrew the septum cap of the injector.
  - b. Unscrew the ring nut.
  - c. Remove the septum holder/liner cap with septum from the injector body head.
8. Replace the head body o-rings.
  - a. Use tweezers to remove the body head internal and external o-rings, and replace them with new o-rings.
9. Reinstall the top parts of the injector.
  - a. Place the septum holder/liner cap with the septum on the body head of the injector and fix it by screwing the ring nut.
  - b. Screw and tighten the injector cap to finger-tight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.


10. Close the module flap cover.
11. If present, move the autosampler towards the module to restore the original alignment.
12. Wait two or three minutes for helium to purge the inlet.
13. Turn the inlet flow back on and enable the inlet heater.
14. Set the normal injector, detector, and GC working conditions.

## Maintaining a Programmable Temperature Vaporizing Injector (PTV)

This section provides instructions for maintaining a Programmable Temperature Vaporizing injector (PTV).

The module and injector components are shown in [Figure 102](#) and [Figure 103](#).


**Figure 102.** PTV Module Components

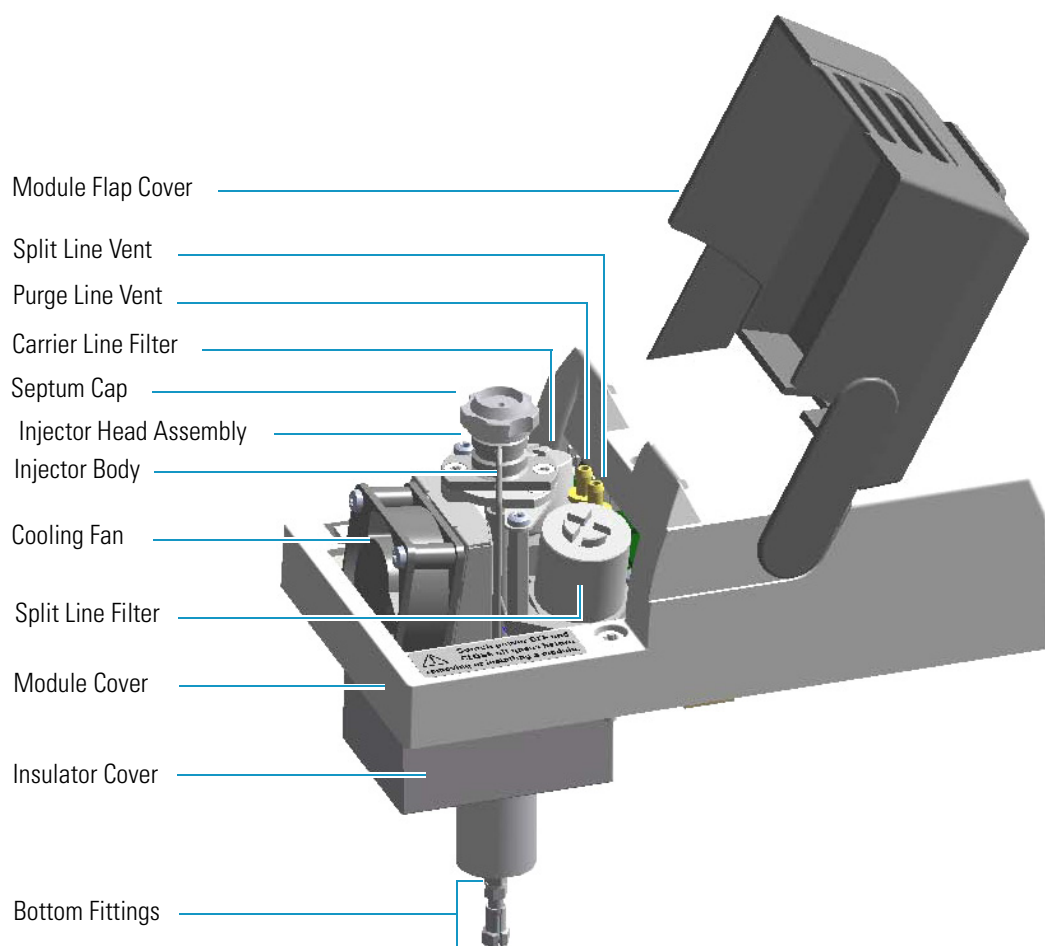

### 3 Performing Injectors Routine Maintenance

Maintaining a Programmable Temperature Vaporizing Injector (PTV)

**Figure 103.** PTV Injector Components

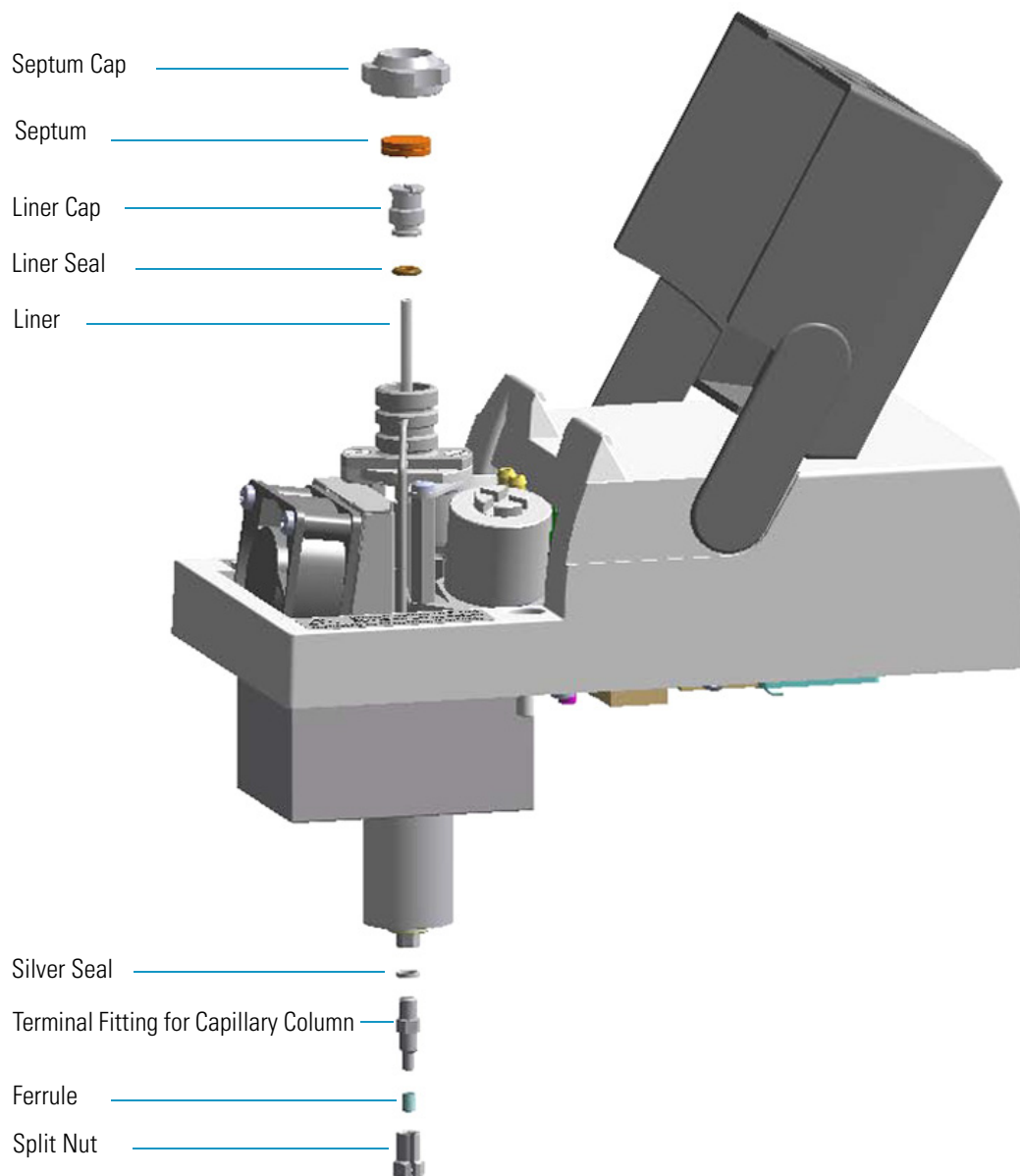

The PTV injector periodic maintenance includes:

- **Replacing the PTV septum**

Replace the septum at least after every 200 injections, or every time a problem related to septum damage, or wear occurs.

See [“Replacing the PTV Septum”](#) on [page 120](#).

- **Cleaning or replacing the PTV liner**

Check the liner for contaminants, debris, breakage and proper installation.

The liner must be replaced depending on the number of injections performed, and the characteristics of the samples injected. Typical symptoms will indicate that the liner must be replaced. The most common is the appearance of tailing peaks in the chromatogram, particularly for polar compounds.

**Tip** It is good practice to replace the septum every time you replace the gas liner.

See “Cleaning or Replacing the PTV Glass Liner” on page 121.

When replacing or removing a glass liner, it might break inside the injector. In this case the broken parts of the liner must be removed from the injector, including the glass splinters that might fall into the lower part of the vaporization chamber.

See “Replacing the PTV Broken Liner” on page 123.

- **Replacing the active carbon filters on the carrier gas line and split line**

The active carbon filters must be replaced depending on the volume of solvent injected in the time.

See “Replacing the PTV Carrier and Split Lines Filters” on page 125.

Before maintaining the injector, read the following warning:



**WARNING** The injector fittings could be hot. Carry out all the operation at low temperature to avoid burns. Therefore, before beginning the sequence, the injector must be cooled to room temperature.



**CAUTION** When handling organic solvents you must take precautions to avoid health hazard.

#### **Materials needed to maintain a PTV injector**

Septum

Tweezers

Glass liner

Liner seal (O-ring)

Ultrasonic cleaner

Mixture 1:1 methanol/acetone

Carrier gas line and/or split gas line active carbon filters

5.5 x 25 mm slotted stubby driver

Silver seal

### 3 Performing Injectors Routine Maintenance

Maintaining a Programmable Temperature Vaporizing Injector (PTV)


**Note** For maintaining or replacing any other component not listed in this section, see Chapter 6, “Injectors Advanced Maintenance.”

- “Baking-out Contaminants from SSL, SSLBKF, HeS-S/SL, PTV, and PTVBKF Injectors” on page 230
- “Removing/Replacing an Injector Module” on page 231
- “Cleaning the PTV Injector Head Assembly” on page 246

## Replacing the PTV Septum

### ❖ To replace the septum

**Figure 104.** PTV Injector: Septum Replacement

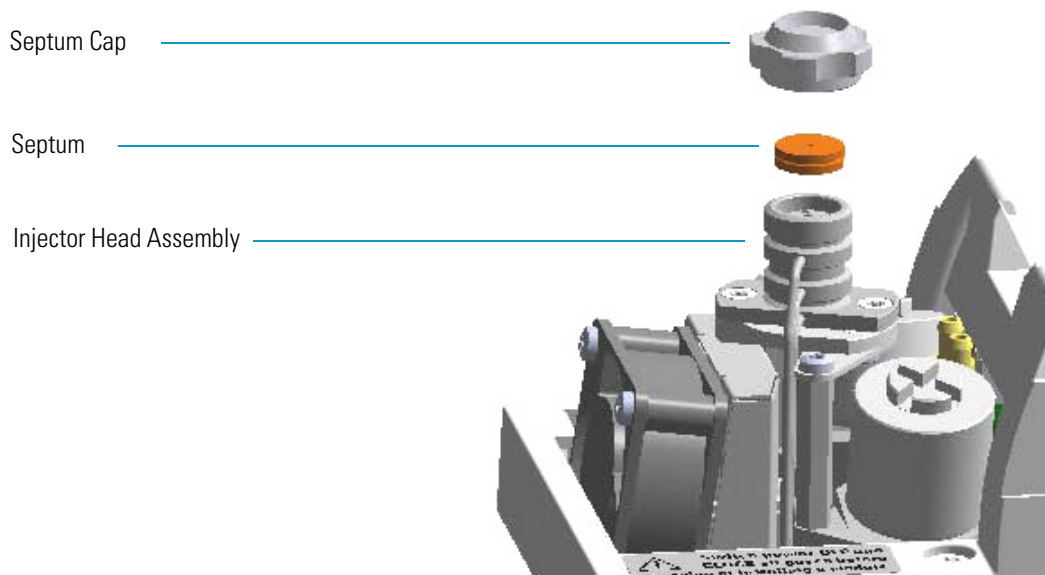

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Replace the septum.
  - a. Unscrew the septum cap of the injector.
  - b. Using tweezers, remove the septum from the injector head assembly.

- c. Avoid touching the septum with your fingers. Insert a new septum into the injector head assembly using tweezers. Screw and tighten the septum cap to finger-tight.




**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

7. Close the module flap cover.
8. If present, move the autosampler towards the module to restore the original alignment.
9. Turn the carrier gas on.
10. Set the normal injector, detector, and GC working conditions.

## Cleaning or Replacing the PTV Glass Liner

### ❖ To clean or replace the glass liner

**Figure 105.** PTV Injector: Glass Liner Replacement

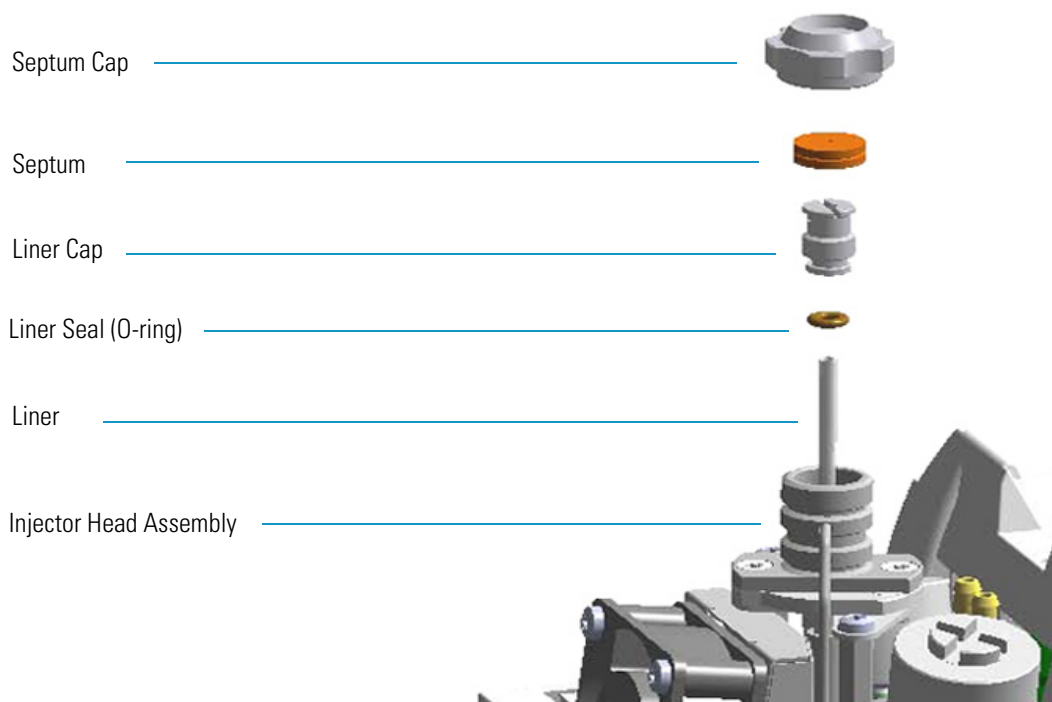

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

### 3 Performing Injectors Routine Maintenance

Maintaining a Programmable Temperature Vaporizing Injector (PTV)

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the top part of the injector.
  - a. Unscrew the septum cap of the injector.
  - b. Using tweezers, remove the septum from the top of the injector head assembly.
7. Remove the liner.
  - a. Using the slotted stubby driver provided, unscrew and remove the liner cap.
  - b. Use tweezers to remove the liner with the liner seal (o-ring) from the injector.



**CAUTION** Be careful not to break the glass liner when removing it. Glass splinters might fall into the lower part of the vaporization chamber. If the glass liner breaks, see [“To replace a broken liner”](#) on [page 86](#).

8. Replace or clean the liner.
  - If you are going to use a new liner, go directly to [step 10](#).
  - If you are going to clean the dirty liner, go to [step 9](#).
9. Clean the liner.
  - a. Put the liner into an ultrasonic cleaner filled with a methanol/acetone mixture (1:1).
  - b. Sonicate the liner for about half an hour.
  - c. Using tweezers, remove the liner from the bath, and dry it with compressed clean air.
10. Install the liner
  - a. Holding the new (or cleaned) liner with tweezers place a new liner seal over the liner.
  - b. Insert the liner into the injector, and push it gently towards the bottom of the injector.
  - c. Reinsert and screw the liner cap using the slotted stubby driver.
11. Reinstall the top parts of the injector.
  - a. Place the septum into the injector head assembly.



**CAUTION** It is suggested to replace the septum with a new one.

- b. Screw and tighten the septum cap to finger-tight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

12. Close the module flap cover.




13. If present, move the autosampler towards the module to restore the original alignment.
14. Turn the carrier gas on.
15. Set the normal injector, detector, and GC working conditions.

## Replacing the PTV Broken Liner

### ❖ To replace a broken liner

**Figure 106.** PTV Injector: Broken Liner Replacement

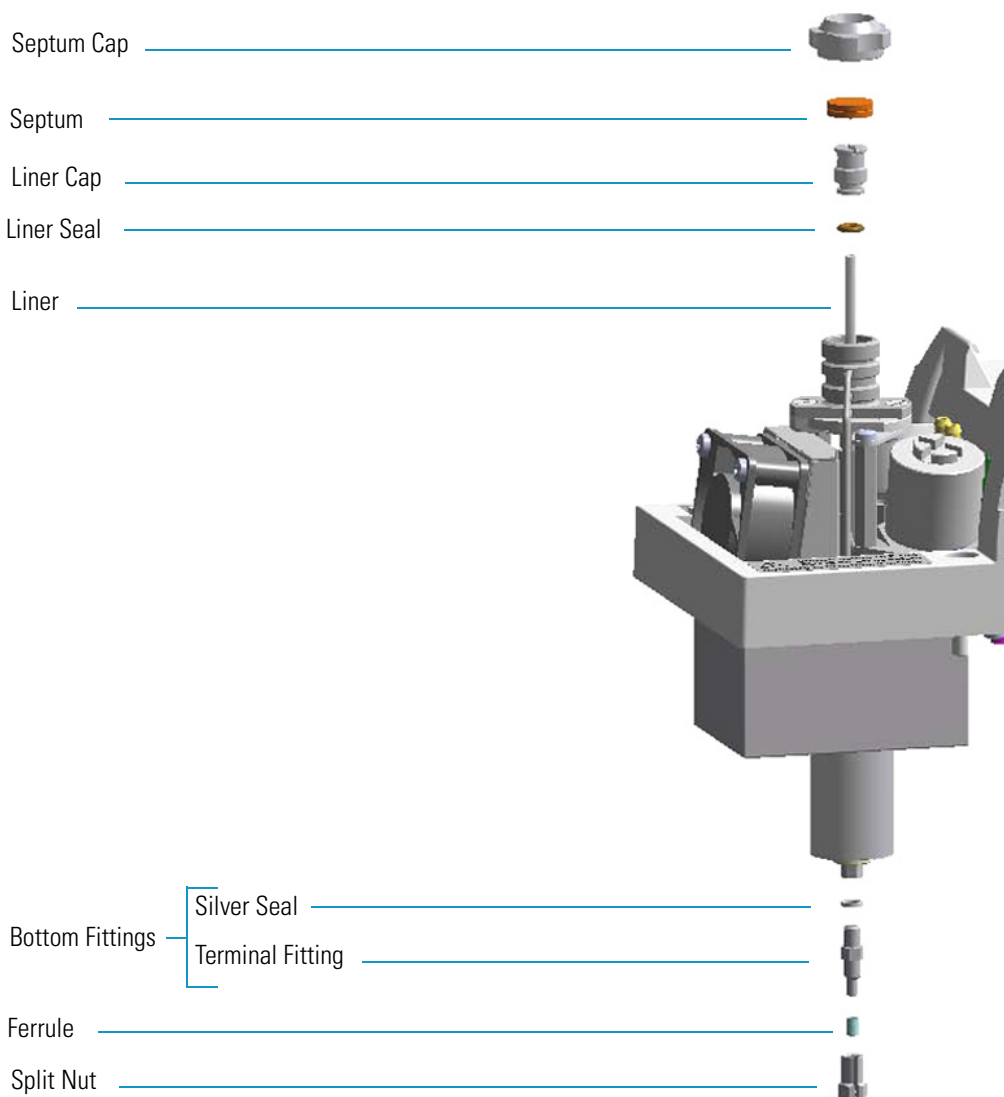

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

### 3 Performing Injectors Routine Maintenance

Maintaining a Programmable Temperature Vaporizing Injector (PTV)

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the top part of the injector.
  - a. Unscrew the septum cap of the injector.
  - b. Using tweezers, remove the septum from the top of the injector head assembly.
7. Remove the broken liner.
  - a. Using the slotted stubby driver provided, unscrew and remove the liner cap.
  - b. Use tweezers to remove the upper part of broken liner with the liner seal from the injector.
8. Remove the bottom parts of the injector.
  - a. Inside the oven, unscrew the split nut, and remove the analytical column with its ferrule from the terminal fitting at the bottom of the injector.
  - b. Unscrew the terminal fitting with the silver seal. Glass splinters from the broken liner will fall from the injector.
  - c. With the aid of a pipe cleaner, remove the possible glass fragments from the vaporization chamber.
9. Reinstall the bottom part of the injector.
  - a. Reinstall the terminal fitting with the silver seal.
  - b. Reinstall the analytical column.
10. Install a new liner.
  - a. Holding the new liner with tweezers, place a new liner seal over the liner.
  - b. Insert the liner into the injector, and push it gently towards the bottom of the injector.
  - c. Reinsert and screw the liner cap using the slotted stubby driver.
11. Reinstall the top parts of the injector.
  - a. Place the septum into the injector head assembly.



**CAUTION** It is suggested to replace the septum with a new one.

- b. Screw and tighten the injector cap to finger-tight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

12. Close the module flap cover.
13. If present, move the autosampler towards the module to restore the original alignment.
14. Turn the carrier gas on.
15. Set the normal injector, detector, and GC working conditions.


## Replacing the PTV Carrier and Split Lines Filters



**IMPORTANT** The dimensions of the filters are different. The filter on the split gas line is bigger than the filter on the carrier gas line. Do not invert their position when you replace them. It is not necessary to replace the filters together.

❖ **To replace the active carbon filters on carrier gas line and split line**

**Figure 107.** PTV Injector: Active Filters Replacement

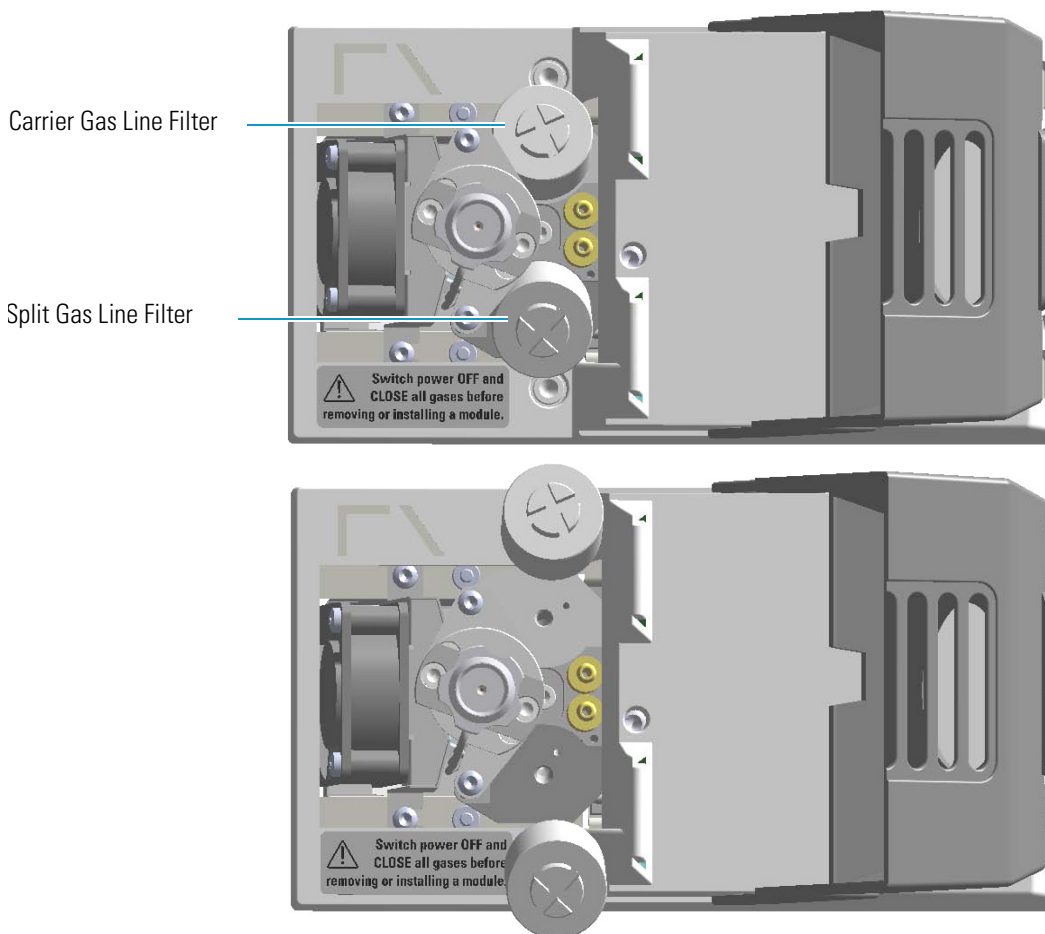

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


4. Put the autosampler away if present.
5. Open the module flap cover.
6. Replace the filter.
  - a. Remove the filter to replace from its seat by turning it counter-clockwise.
  - b. Install the new filter, with the o-rings, in its seat by turning it clockwise.
7. Close the module flap cover.
8. If present, move the autosampler towards the module to restore the original alignment.
9. Turn the carrier gas on.
10. Set the normal injector, detector, and GC working conditions.

## Maintaining a Programmable Temperature Vaporizing Injector with Backflush (PTVBKF)

This section provides instructions for maintaining a Programmable Temperature Vaporizing injector for backflush (PTVBKF) applications.

The module and injector components are shown in [Figure 108](#) and [Figure 109](#).


**Figure 108.** PTVBKF Module Components

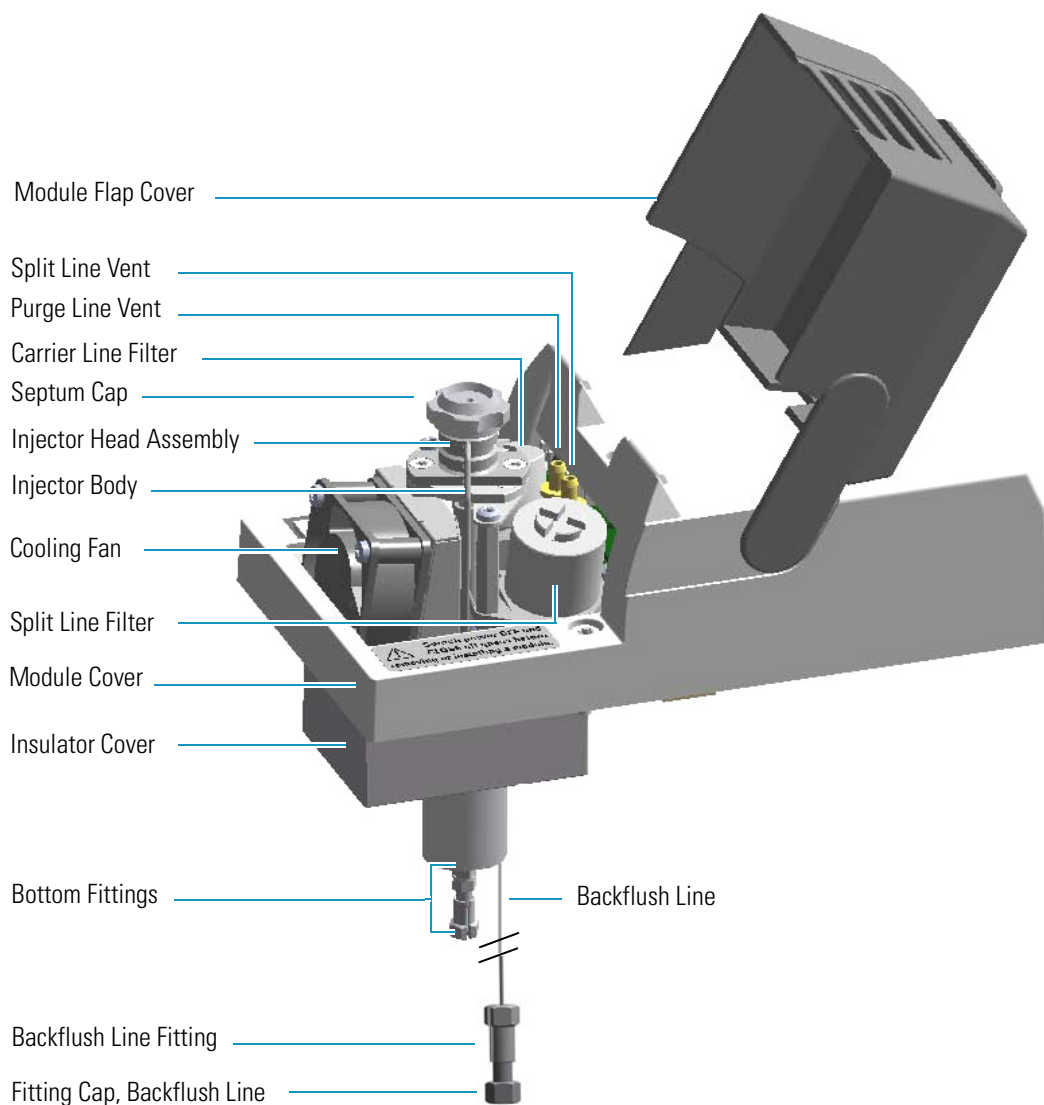

### 3 Performing Injectors Routine Maintenance

Maintaining a Programmable Temperature Vaporizing Injector with Backflush (PTVBKF)

**Figure 109.** PTVBKF Injector Components

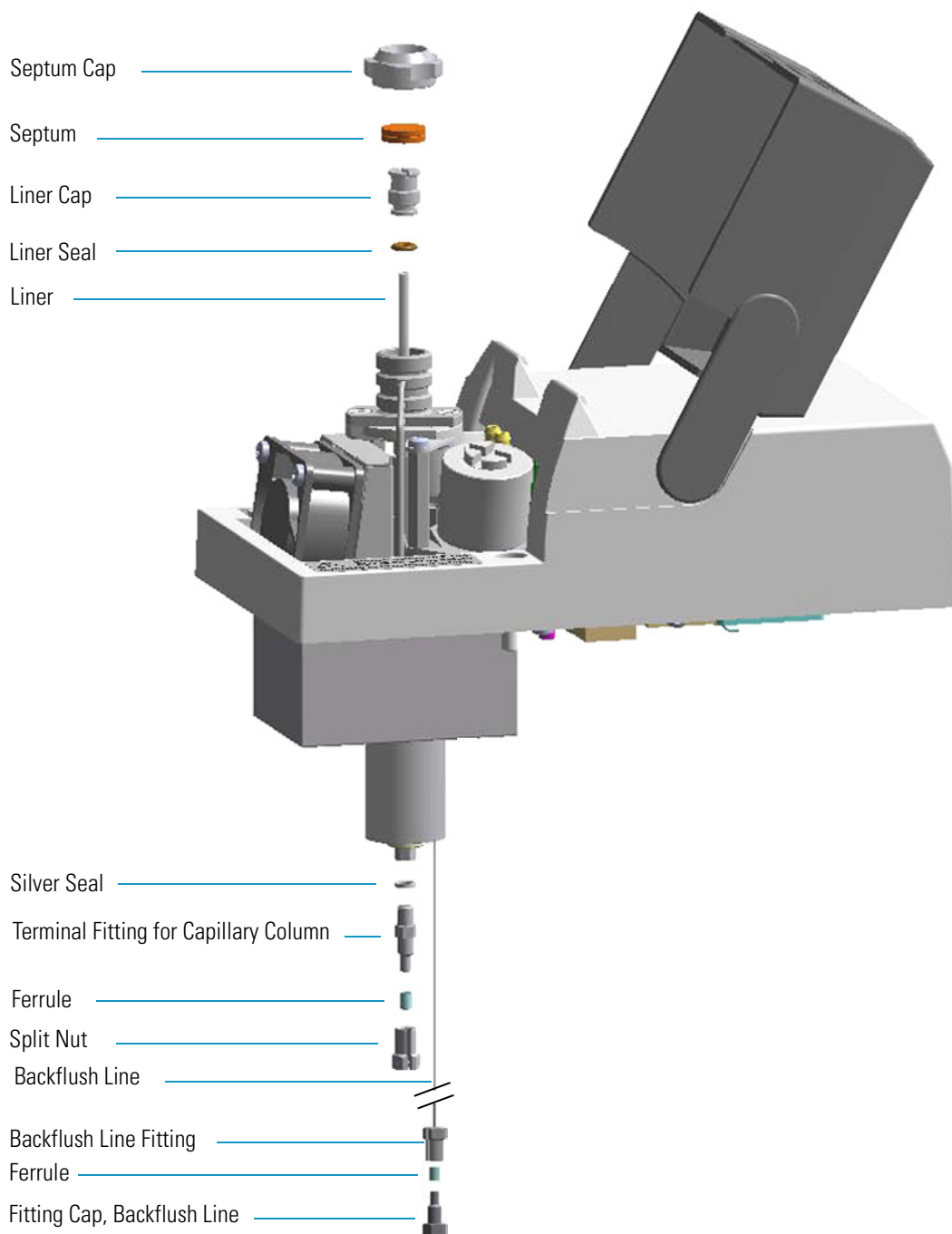

The PTVBKF injector periodic maintenance includes:

- **Replacing the PTVBKF septum**

Replace the septum at least after every 200 injections, or every time a problem related to septum damage, or wear occurs.

See [“Replacing the PTVBKF Septum”](#) on [page 130](#).

- **Cleaning or replacing the PTVBKF liner**

Check the liner for contaminants, debris, breakage and proper installation. The liner must be replaced depending on the number of injections performed, and the characteristics of the samples injected. Typical symptoms will indicate that the liner must be replaced. The most common is the appearance of tailing peaks in the chromatogram, particularly for polar compounds.

See [“Cleaning or Replacing the PTVBKF Glass Liner”](#) on [page 131](#)

**Tip** It is good practice to replace the septum every time the liner is replaced.

When replacing or removing a glass liner, it might break inside the injector. In this case the broken parts of the liner must be removed from the injector, including the glass splinters that might fall into the lower part of the vaporization chamber.

See [“Replacing the PTVBKF Broken Liner”](#) on [page 133](#).

- **Replacing the PTVBKF active carbon filters on the carrier gas line and split line**

The active carbon filters must be replaced depending on the volume of solvent injected in the time.

See [Replacing the PTVBKF Carrier and Split Lines Filters](#) on [page 135](#).

Before maintaining the injector, read the following warning:



**WARNING** The injector fittings could be hot. Carry out all operations at low temperature to avoid burns. Therefore, before beginning the sequence, the injector must be cooled to room temperature.



**CAUTION** When handling organic solvents you must take precautions to avoid health hazards.

#### Materials needed to maintain a PTVBKF injector

Septum

Tweezers

Glass liner

Liner seal (O-ring)

Ultrasonic cleaner

Mixture 1:1 methanol/acetone

Carrier gas line and/or split gas line active carbon filters

5.5 x25 mm slotted stubby driver

Silver seal

### 3 Performing Injectors Routine Maintenance

Maintaining a Programmable Temperature Vaporizing Injector with Backflush (PTVBKF)


**Note** For maintaining or replacing any other component not listed in this section, see Chapter 6, “Injectors Advanced Maintenance.”

- “Baking-out Contaminants from SSL, SSLBKF, HeS-S/SL, PTV, and PTVBKF Injectors” on page 230
- “Removing/Replacing an Injector Module” on page 231
- “Cleaning the PTVBKF Injector Head Assembly” on page 253

## Replacing the PTVBKF Septum

### ❖ To replace the septum

**Figure 110.** PTVBKF Injector: Septum Replacement

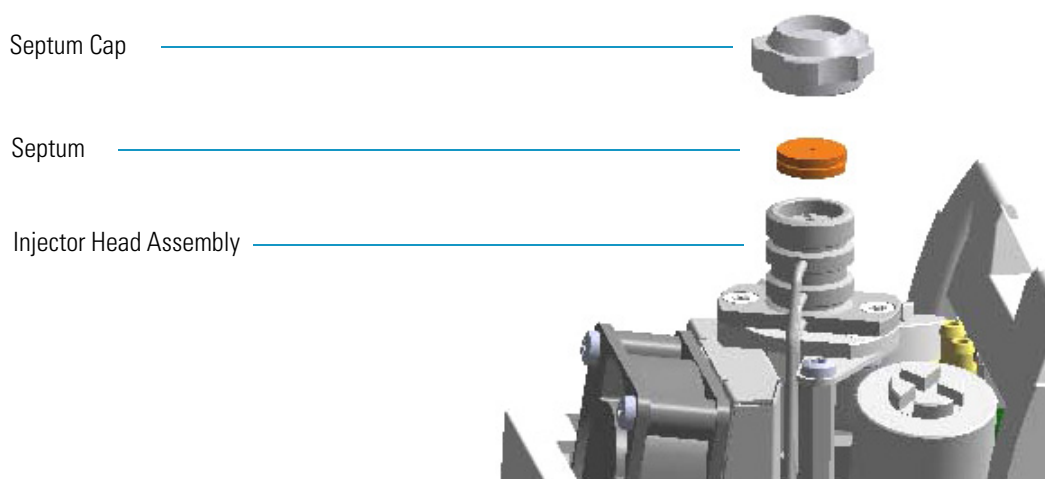

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Replace the septum.
  - a. Unscrew the septum cap of the injector.
  - b. Using tweezers, remove the septum from the injector head assembly.
  - c. Avoid touching the septum with your fingers. Insert a new septum into the injector head assembly using tweezers. Screw and tighten the septum cap to finger-tight.






**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

7. Close the module flap cover.
8. If present, move the autosampler towards the module to restore the original alignment.
9. Turn the carrier gas on.
10. Set the normal injector, detector, and GC working conditions.

## Cleaning or Replacing the PTVBKF Glass Liner

### ❖ To clean or replace the glass liner

**Figure 111.** PTVBKF Injector: Glass Liner Replacement

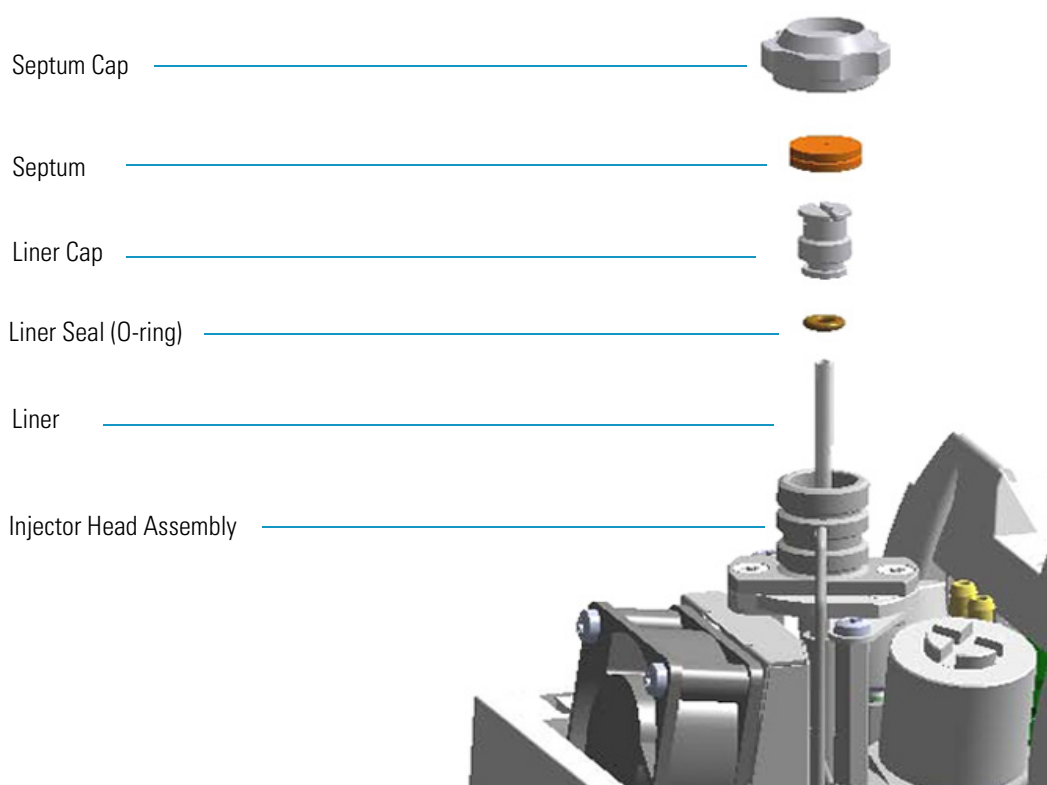

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

### 3 Performing Injectors Routine Maintenance

Maintaining a Programmable Temperature Vaporizing Injector with Backflush (PTVBKF)

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the top part of the injector.
  - a. Unscrew the septum cap of the injector.
  - b. Using tweezers, remove the septum from the top of the injector head assembly.
7. Remove the liner.
  - a. Using the slotted stubby driver provided, unscrew and remove the liner cap.
  - b. Use tweezers to remove the liner with the liner seal (o-ring) from the injector.



**CAUTION** Be careful not to break the glass liner when removing it. Glass splinters might fall into the lower part of the vaporization chamber. If the glass liner breaks, see “Replacing the PTVBKF Broken Liner” on page 133.

8. Replace or clean the liner.
  - If you are going to use a new liner, go directly to [step 9](#).
  - If you are going to clean the dirty liner, go to [step 10](#).
9. Clean the liner.
  - a. Put the liner into an ultrasonic cleaner filled with a methanol/acetone mixture (1:1).
  - b. Sonicate the liner for about half an hour.
  - c. Using tweezers, remove the liner from the bath, and dry it with compressed clean air.
10. Install the liner.
  - a. Holding the new (or cleaned) liner with tweezers, place a new liner seal over the liner.
  - b. Insert the liner into the injector, and push it gently towards the bottom of the injector.
  - c. Reinsert and screw the liner cap using the slotted stubby driver.
11. Reinstall the top parts of the injector.
  - a. Place the septum into the injector head assembly.



**CAUTION** It is suggested to replace the septum with a new one.

- b. Screw and tighten the septum cap to finger-tight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.


12. Close the module flap cover.

13. If present, move the autosampler towards the module to restore the original alignment.
14. Turn the carrier gas on.
15. Set the normal injector, detector, and GC working conditions.

## Replacing the PTVBKF Broken Liner

### ❖ To replace a broken liner

**Figure 112.** PTVBKF Injector: Broken Liner Replacement

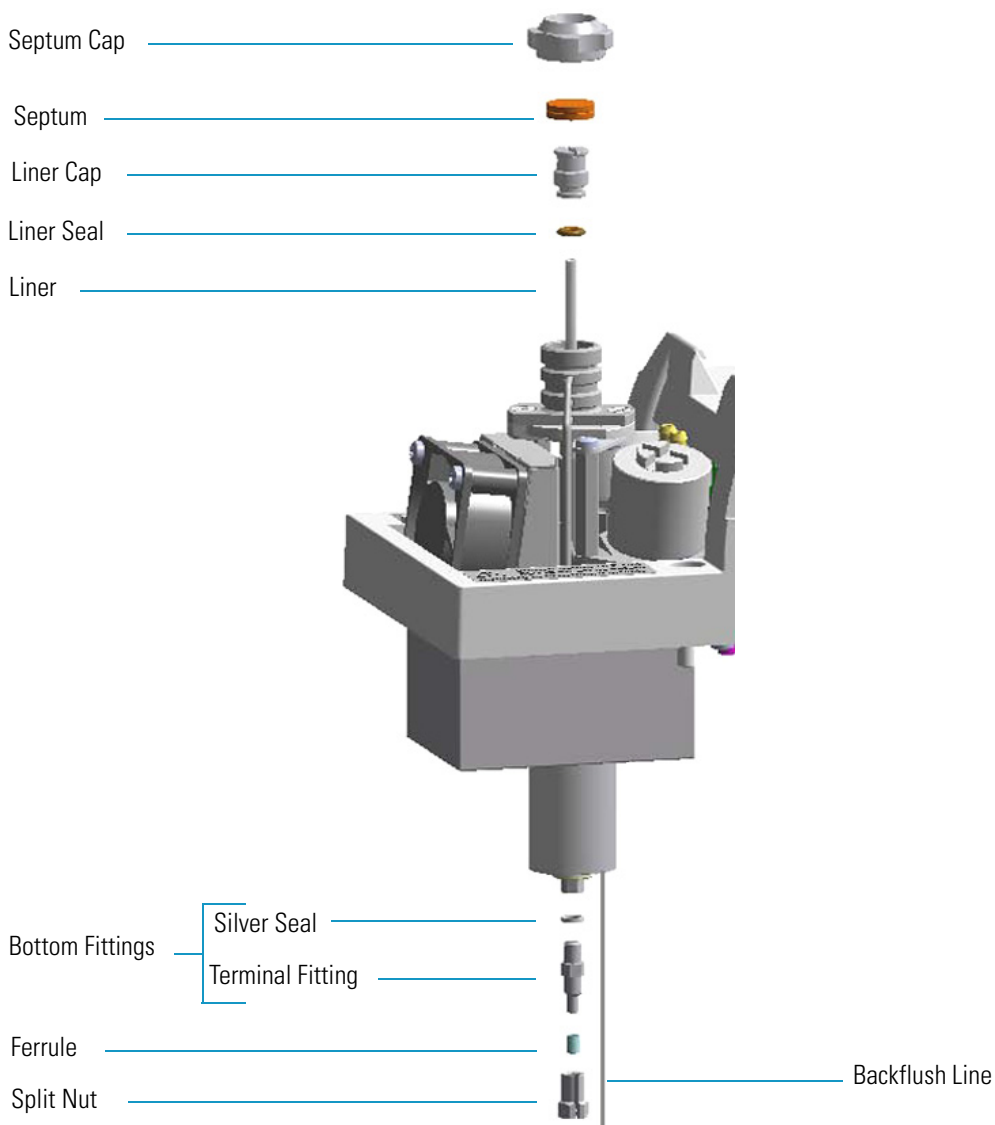

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

### 3 Performing Injectors Routine Maintenance

Maintaining a Programmable Temperature Vaporizing Injector with Backflush (PTVBKF)

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the top part of the injector.
  - a. Unscrew the septum cap of the injector.
  - b. Using tweezers, remove the septum from the top of the injector head assembly.
7. Remove the broken liner.
  - a. Using the slotted stubby driver provided, unscrew and remove the liner cap.
  - b. Use tweezers to remove the upper part of broken liner with the liner seal from the injector.
8. Remove the bottom parts of the injector.
  - a. Inside the oven, unscrew the split nut, and remove the analytical column with its ferrule from the terminal fitting at the bottom of the injector.
  - b. Unscrew the terminal fitting with the silver seal. Glass splinters from the broken liner will fall from the injector.
  - c. With the aid of a pipe cleaner, remove the possible glass fragments from the vaporization chamber.
9. Reinstall the bottom part of the injector.
  - a. Reinstall the terminal fitting with the silver seal (if necessary replace it with a new one).
  - b. Reinstall the analytical column.
10. Install a new liner.
  - a. Holding the new liner with tweezers, place a new liner seal over the liner.
  - b. Insert the liner into the injector, and push it gently towards the bottom of the injector.
  - c. Reinsert and screw the liner cap using the slotted stubby driver.
11. Reinstall the top parts of the injector.
  - a. Place the septum into the injector head assembly.



**CAUTION** It is suggested to replace the septum with a new one.

- b. Screw and tighten the injector cap to finger-tight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.


12. Close the module flap cover.
13. If present, move the autosampler towards the module to restore the original alignment.
14. Turn the carrier gas on.
15. Set the normal injector, detector, and GC working conditions.

## Replacing the PTVBKF Carrier and Split Lines Filters

**IMPORTANT** The dimensions of the filters are different. The filter on the split gas line is bigger than the filter on the carrier gas line. Do not invert their position when you replace them. It is not necessary to replace the filters together.

### ❖ To replace the active carbon filters on carrier gas line and split line

**Figure 113.** PTVBKF Injector: Active Filters Replacement

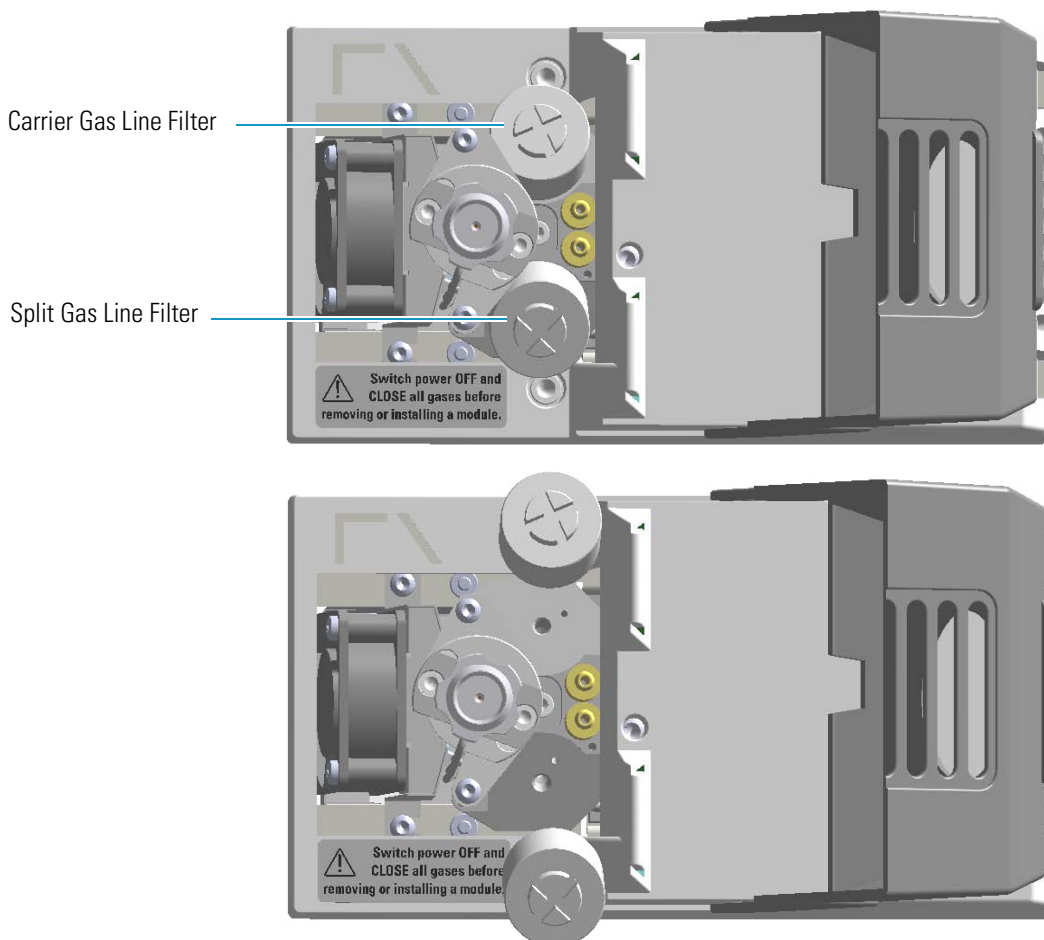

### 3 Performing Injectors Routine Maintenance

Maintaining a Programmable Temperature Vaporizing Injector with Backflush (PTVBKF)

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Replace the filter.
  - a. Remove the filter to replace from its seat by turning it counter-clockwise.
  - b. Install the new filter, with the o-rings, in its seat by turning it clockwise.
7. Close the module flap cover.
8. If present, move the autosampler towards the module to restore the original alignment.
9. Turn the carrier gas on.
10. Set the normal injector, detector, and GC working conditions.

## Performing Detectors Routine Maintenance

This chapter provides instructions for performing routine maintenance on the TRACE 1300 and TRACE 1310 detectors modules.


### Contents

- [Maintaining a Flame Ionization Detector \(FID\)](#)
- [Maintaining a Nitrogen Phosphorous Detector \(NPD\)](#)
- [Maintaining a Thermal Conductivity Detector \(TCD\)](#)
- [Maintaining an Electron Capture Detector \(ECD\)](#)
- [Maintaining a Flame Photometric Detector \(FPD\)](#)
- [Maintaining a Pulsed Discharge Detector \(PDD\)](#)

## Maintaining a Flame Ionization Detector (FID)


This section provides instructions for maintaining a Flame Ionization Detector (FID). The module and detector components are shown in [Figure 114](#) and [Figure 115](#).

**Figure 114.** FID: Module Components

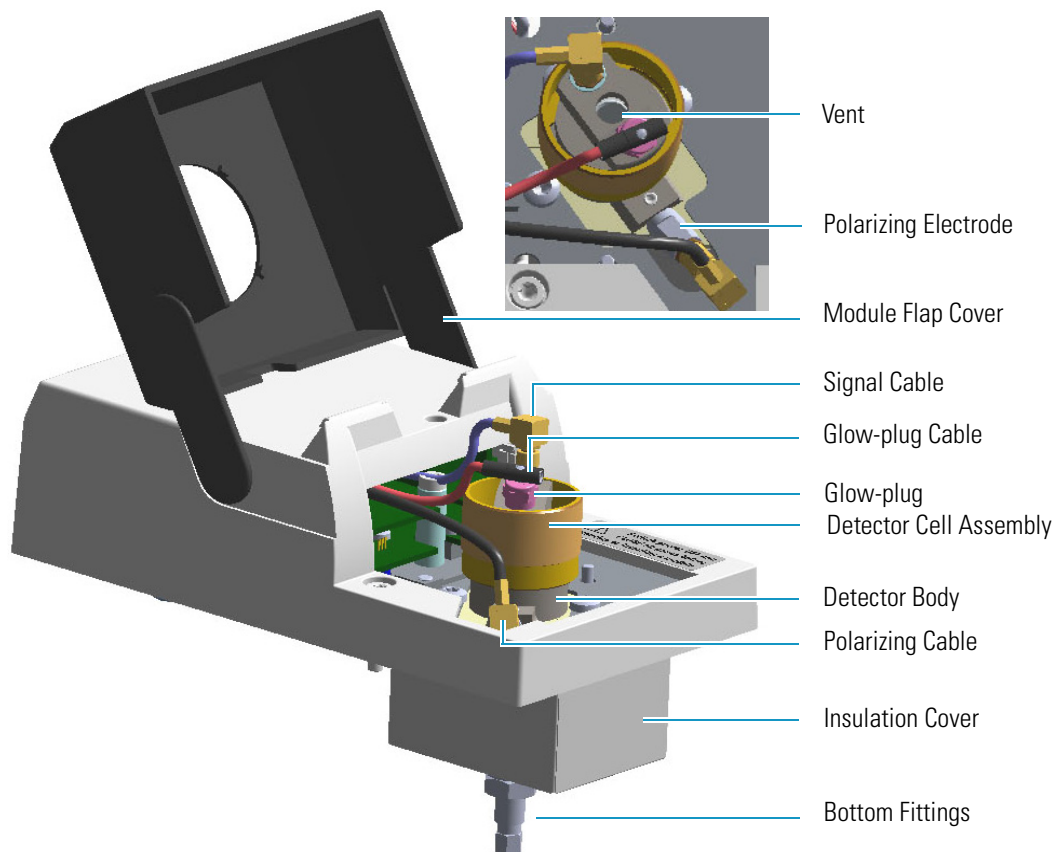




**Figure 115.** FID: Cell Assembly Components

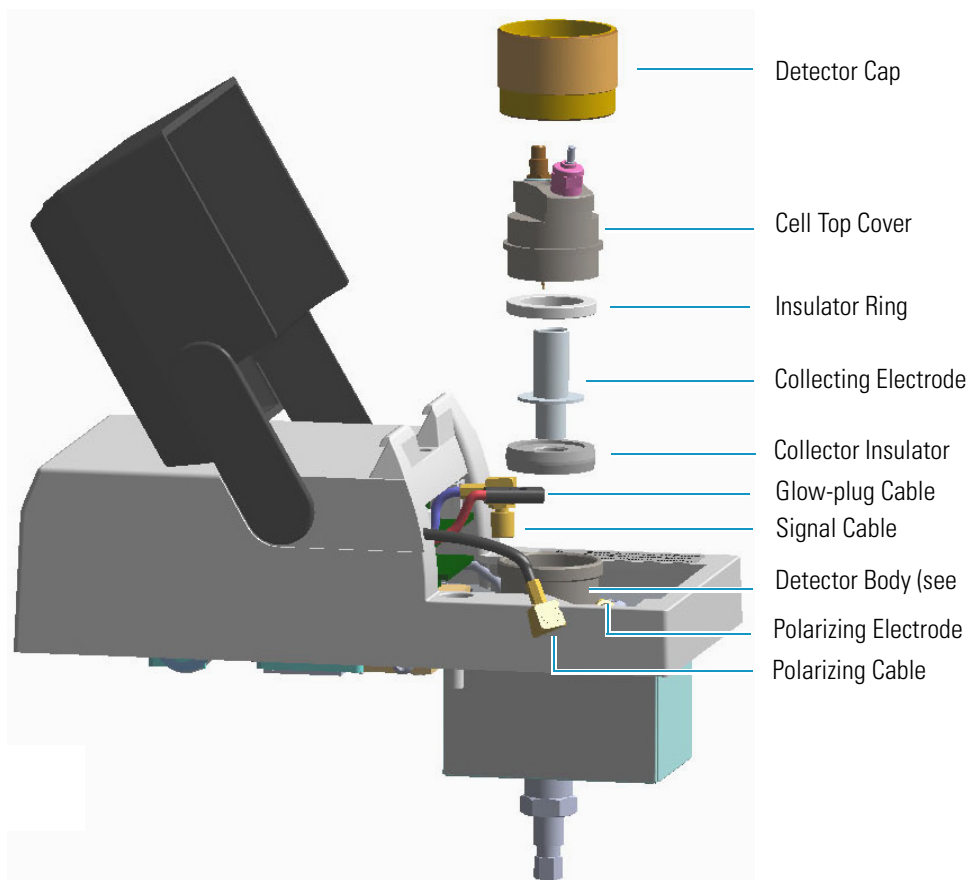

**Figure 116.** FID: Jet and Polarizing Electrode

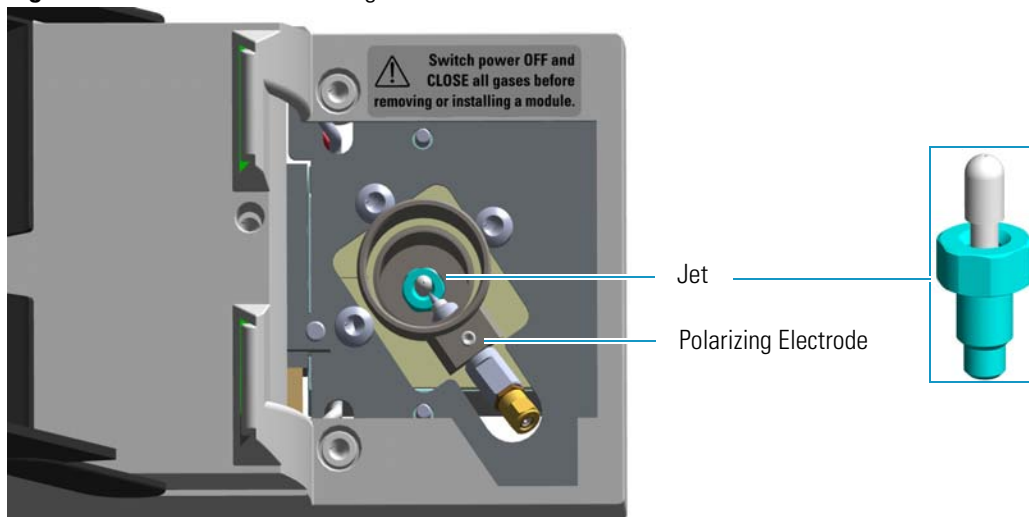

To keep optimum performance of the FID, you must keep it clean and free of dust and deposits. Symptoms, such as reduced sensitivity and increased noise, indicate that the detector needs cleaning.

## 4 Performing Detectors Routine Maintenance

### Maintaining a Flame Ionization Detector (FID)

The FID detector periodic maintenance includes:

- **Cleaning the collecting electrode**

It is suggest to clean the collecting electrode annually, or when a high noise baseline is found, due to some compounds that produce detector contamination. If necessary, replace the collecting electrode.

See “[Cleaning or Replacing the FID Collecting Electrode](#)” on [page 141](#).

- **Replacing the FID ignition glow-plug**

This operation is NOT part on the ordinary maintenance. The flame ignition element must be replaced only when defective.

See “[Replacing the FID Ignition Glow-plug](#)” on [page 144](#).

Before maintaining the detector, read the following precautions:



**WARNING** Carry out all the operations at low temperature to avoid burns. Therefore, before beginning maintenance, cool the detector to room temperature.



**CAUTION** When handling organic solvents you must take precautions to avoid health hazards.

#### Materials needed to maintain a FID detector

Ultrasonic cleaner

Liquid detergent

GC-grade methanol

Distilled water

FID Ignition Glow-plug

Collecting electrode (if necessary)

8-mm elbowed box wrench

T6 Torxhead screwdriver

Forceps or tweezers

**Note** For maintaining or replacing any other component not listed in this section, see [Chapter 7, “Detectors Advanced Maintenance.”](#)

- “[Removing/Replacing a Detector Module](#)” on [page 262](#)
- “[Cleaning or Replacing the FID Jet](#)” on [page 271](#)

## Cleaning or Replacing the FID Collecting Electrode


### ❖ To clean or replace the collecting electrode

1. Put the GC in standby condition.
2. Turn the flame off. The fuel gases, hydrogen and air, are automatically closed.
3. Cool the detector to room temperature.
4. Turn the make-up gas off.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

5. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
6. If external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.
7. Open the module flap cover.
8. Disconnect the signal, glow-plug, and polarizing cables from their contacts on the cell.


**Figure 117.** FID Cables

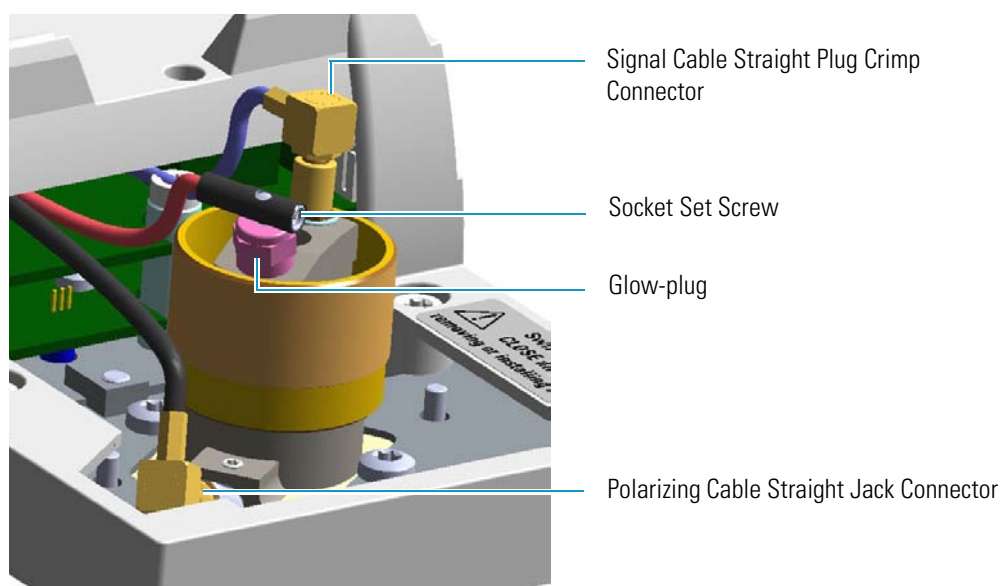

## 4 Performing Detectors Routine Maintenance


### Maintaining a Flame Ionization Detector (FID)

**Figure 118.** FID Cables Disconnection

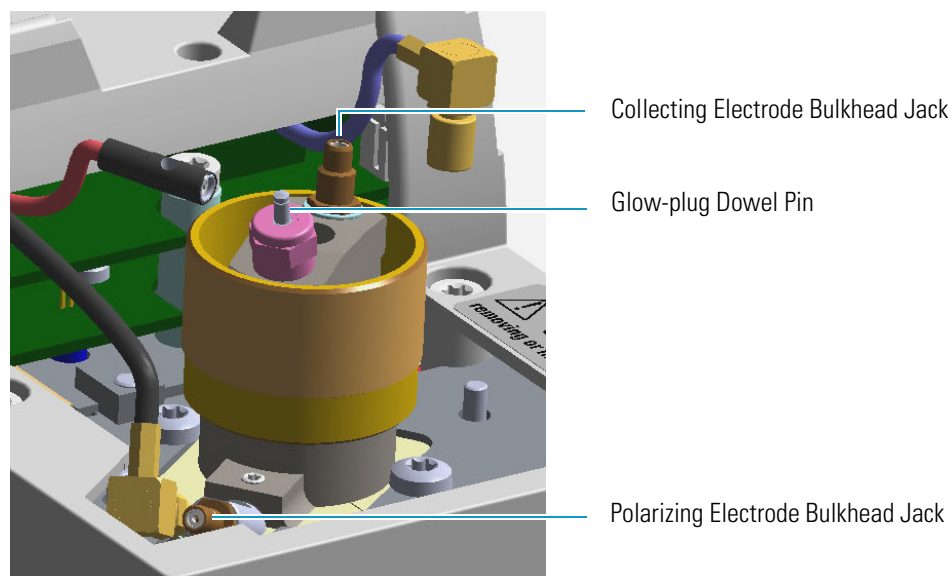

- a. Loosen the glow-plug cable socket set screw using a T6 Torxhead screwdriver. Carefully pull out the terminal body of the glow-plug cable from the dowel pin on the glow-plug element.
  - b. Unscrew and pull out the straight plug crimp connector of the signal cable from the collecting electrode bulkhead jack.
  - c. Pull out the quick coupling straight jack connector of the polarizing cable from the polarizing electrode bulkhead jack.
  - d. Carefully move the cables in order to have free space for handling the detector.
9. Remove the top parts of the detector.

**Figure 119.** FID Collecting Electrode

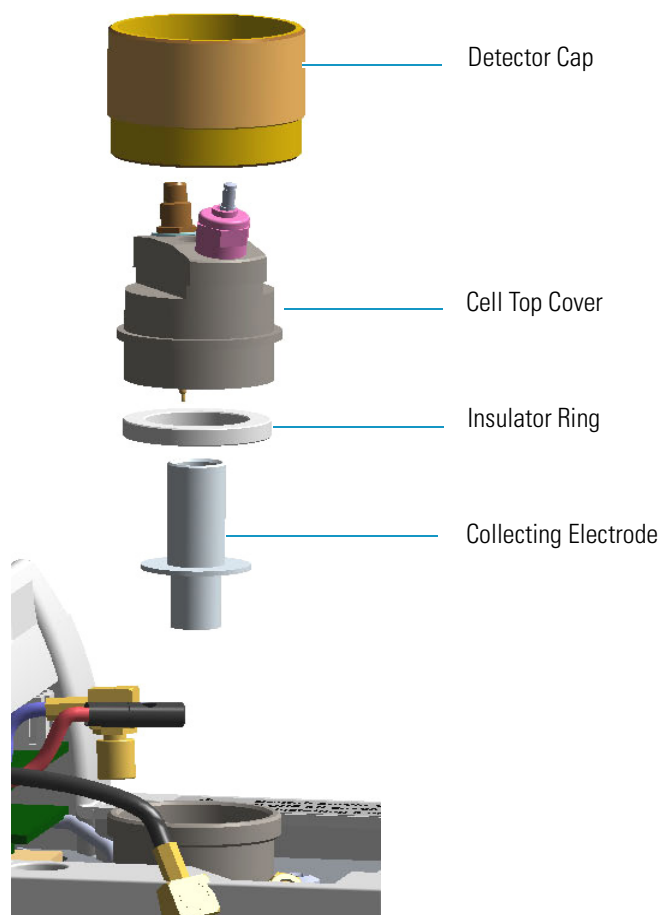

- a. Unscrew and remove the detector cap, paying attention to not rotate the cell top cover.
  - b. Remove the cell top cover and put it on a safe place. Pay attention to not damaging the FID collecting electrode pin.
10. Remove, clean and reinstall the collecting electrode.
- a. Using forceps or tweezers, carefully extract the collecting electrode and the insulator ring through the top of the detector cell.
  - b. Slip off the insulator ring from the collecting electrode.



**WARNING** The insulator ring must be placed carefully away in a place where it stays clean.

- c. Place the collecting electrode in the ultrasonic cleaner filled with liquid detergent.
- d. Sonicate the collecting electrode for about five minutes.
- e. Handling the collecting electrode with forceps or tweezers, rinse it using distilled water, then methanol.

- f. Place the collecting electrode on a paper towel and let it air dry.

**Note** If after the cleaning the physical condition of the collecting electrode does not permit its reuse, replace the electrode with a new one.

- g. Place the cleaned (or a new) collecting electrode and the insulator ring in its housing.
11. Reinstall the top parts of the detector.
  - a. Remount the cell top cover, then screw the detector cap, paying attention to not rotate the cell top cover.
  - b. Reconnect the signal, glow-plug, and polarizing cables.
12. Close the module flap cover.
13. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
14. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
15. Set the make-up gas on.
16. Set the normal detector working conditions.
17. Ignite the flame.

## Replacing the FID Ignition Glow-plug

### ❖ To replace the ignition glow-plug


1. Put the GC in standby condition.
2. Turn the flame off. The fuel gases, hydrogen and air, are automatically closed.
3. Cool the detector to room temperature.
4. Turn the make-up gas off.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


5. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.

6. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
7. Open the module flap cover.
8. Disconnect the signal cable, glow-plug and polarizing cables from their contacts on the cell top cover.

**Figure 120.** FID Cables

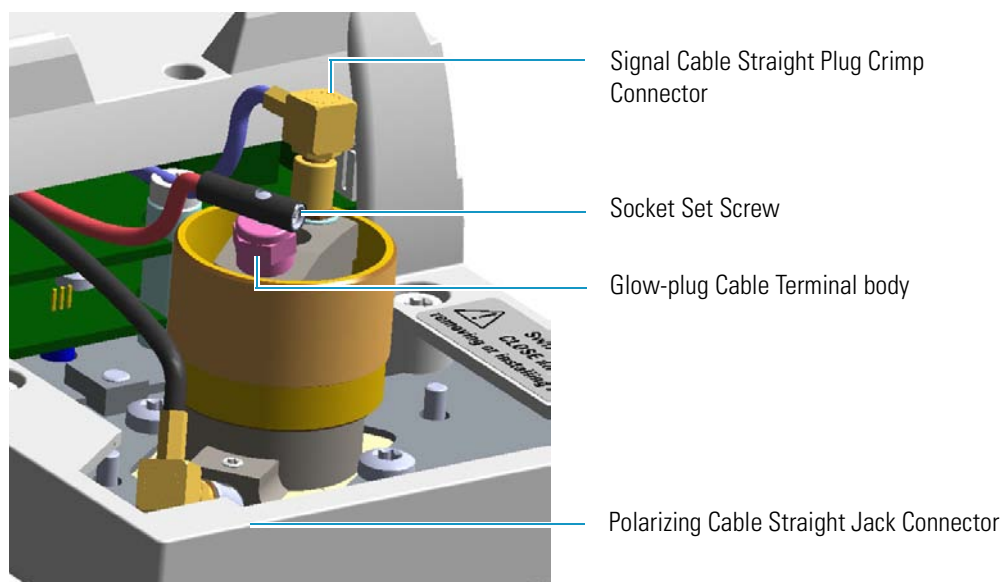

**Figure 121.** FID Cables Disconnection

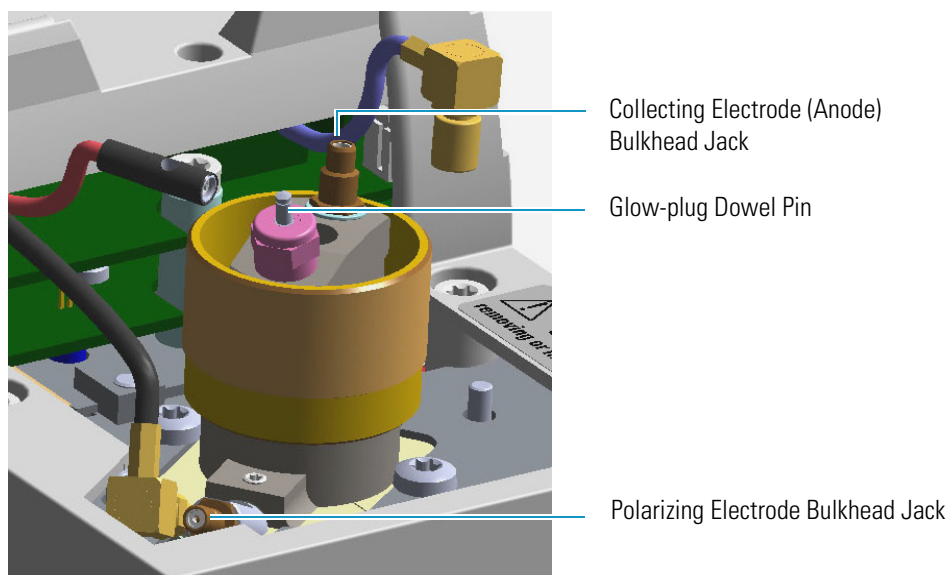


- a. Loosen the glow-plug cable socket set screw using a T6 Torxhead screwdriver. Carefully pull out the terminal body of the glow-plug cable from the dowel pin on the glow-plug element.

## 4 Performing Detectors Routine Maintenance

### Maintaining a Flame Ionization Detector (FID)

- b. Unscrew and pull out the straight plug crimp connector of the signal cable from the collecting electrode (anode) bulkhead jack.
  - c. Pull out the quick coupling straight jack connector of the polarizing cable from the polarizing electrode bulkhead jack.
  - d. Carefully move the cables in order to have free space for handling the detector.
9. Replace the defective glow-plug.

**Figure 122.** FID Glow-plug Replacement

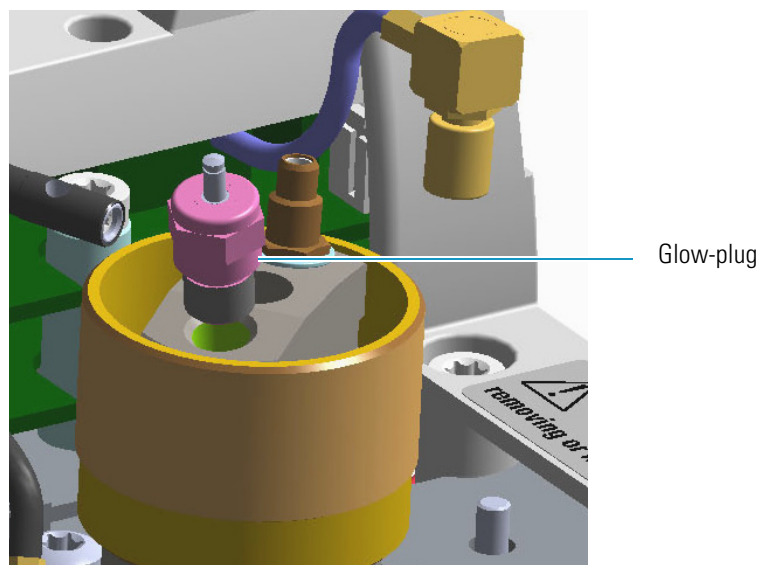


- a. Using an 8-mm wrench, unscrew and remove the defective glow-plug with its washer.
  - b. Replace the glow-plug with a new one interposing the new washer.
  - c. Reconnect the glow-plug, signal, and polarizing cables.
10. Close the module flap cover.
  11. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
  12. Power on the GC.
    - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
    - b. Flip up the power switch (breaker) to the position I.
  13. Turn the make-up gas on.
  14. Set the normal detector working conditions.
  15. Ignite the flame.




## Maintaining a Nitrogen Phosphorous Detector (NPD)

This section provides instructions for maintaining a Nitrogen Phosphorous Detector (NPD). The components of the module and detector are shown in [Figure 123](#), [Figure 124](#), and [Figure 125](#).


**Figure 123.** NPD Module Components

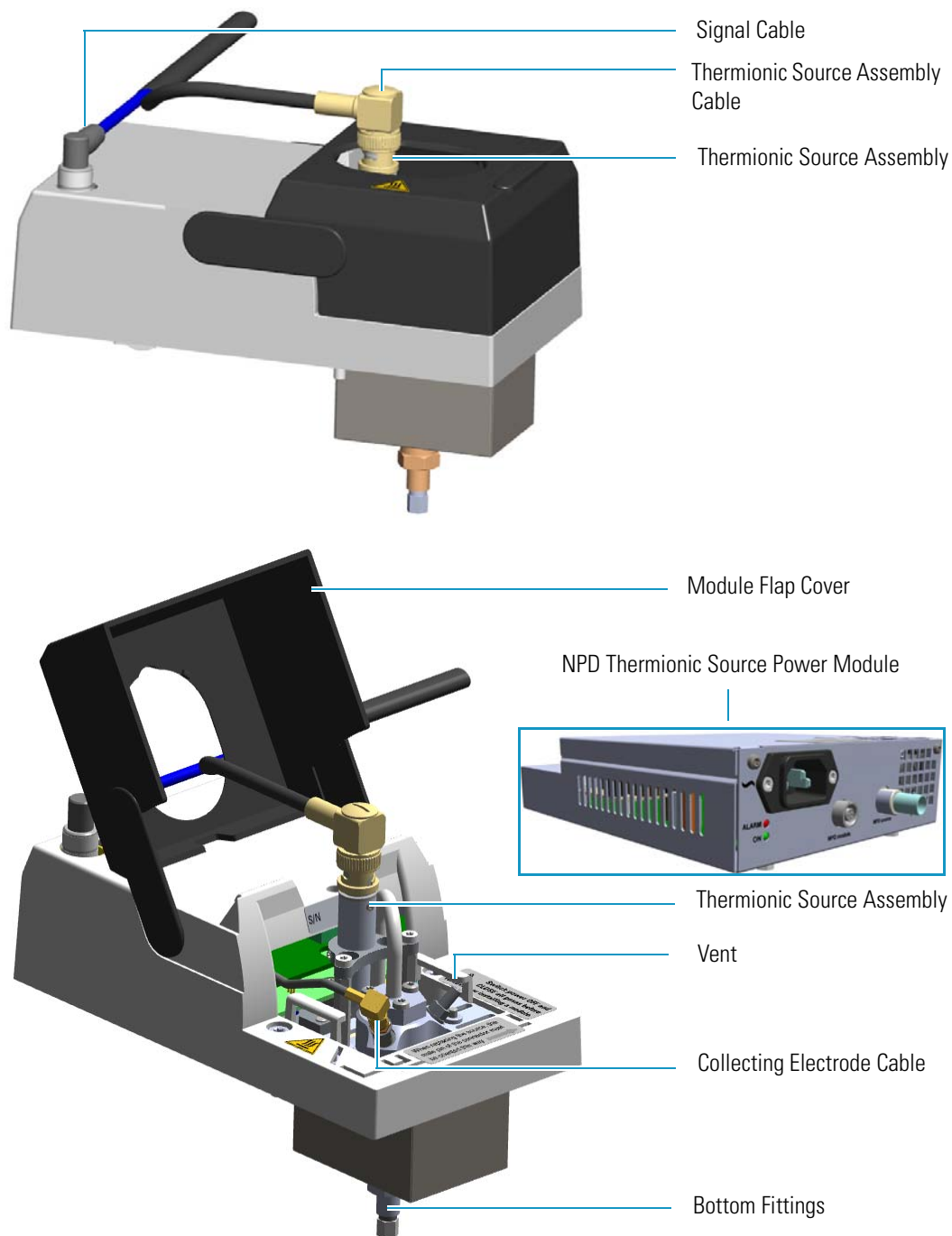

**Figure 124.** NPD Cell Assembly Components

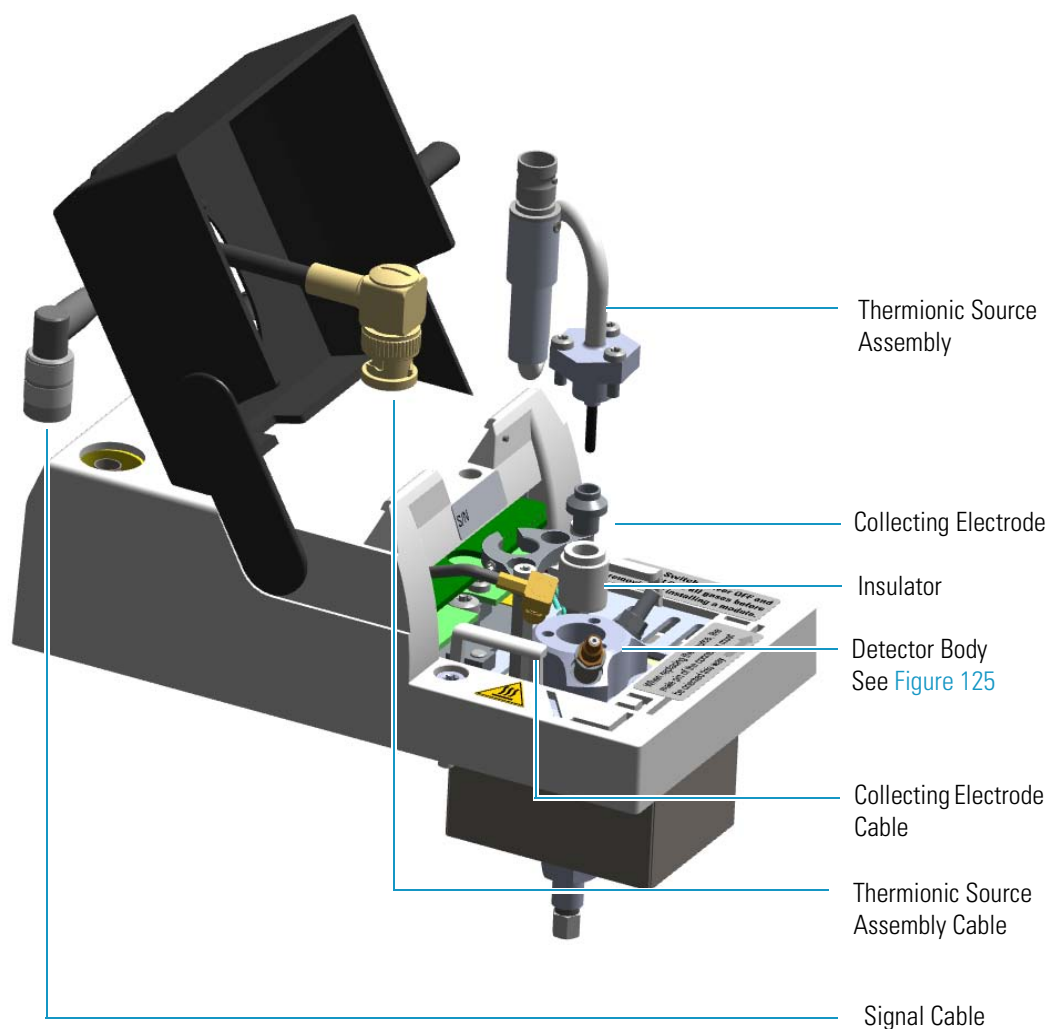

**Figure 125.** NPD Detector: Jet and Collecting Electrode (Anode) Pin

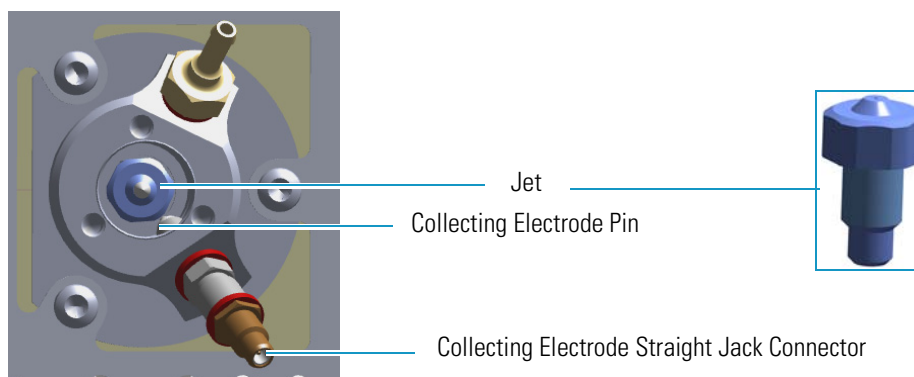

To ensure optimum performance of the NPD, you must keep it clean and free of dust and deposits. Symptoms such as reduced sensitivity and increased noise indicate that detector cleaning, or thermionic source replacement could be necessary.

**Note** The thermionic source is supplied by a separated NPD Thermionic Source Power Module placed into an external module housing on the back of the GC.

The NPD detector periodic maintenance includes:

- **Replacing the thermionic source**

New thermionic sources may require slightly different values of source heating current to produce the same signal observed with a previous thermionic source.

See [“Replacing the NPD Thermionic Source”](#) on [page 150](#).

- **Cleaning the collecting electrode**

We suggest cleaning the collecting electrode semiannually, or when a high noise baseline is found, due to some compounds that produce detector contamination. If necessary, replace the collecting electrode.

See [“Cleaning or Replacing the NPD Collecting Electrode”](#) on [page 156](#).

Before maintaining the detector, read the following precautions:



**WARNING** Carry out all the operations at low temperature to avoid burns. Therefore, before beginning maintenance, cool the detector to room temperature.



**CAUTION** When handling organic solvents you must take precautions to avoid health hazards.

#### **Materials needed to maintain a NPD detector**

Thermionic Source

Gloves

Ultrasonic cleaner

Liquid detergent

GC-grade methanol

Distilled water

Collecting electrode (if necessary)

Forceps or tweezers

T6 Torxhead key

T10 Torxhead screwdriver

6 mm wrench

**Note** For maintaining or replacing any other component not listed in this section, see [Chapter 7, “Detectors Advanced Maintenance.”](#)

- “Measuring the NPD Gas Flows” on page 276
- “Cleaning or Replacing the NPD Jet” on page 285

## Replacing the NPD Thermionic Source



**CAUTION** The thermionic source is delicate. Be careful not to break or crack the source. When performing maintenance on the NPD, avoid touching the source with your fingers, and prevent it from coming in contact with other surfaces. Wear clean, lint-free gloves to prevent contamination of parts with dirt and skin oils.

### ❖ To replace the thermionic source

1. Put the GC in standby condition.
2. Switch off the thermionic source. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.
4. Turn the make-up gas off.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

5. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
6. Unplug the power cable from the AC input connector of the **NPD Thermionic Source Power Module** and from the wall outlet.

If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.


7. Unplug the signal cable from the detector module. See [Figure 126](#).

**Figure 126.** Unplug NPD Signal Cable

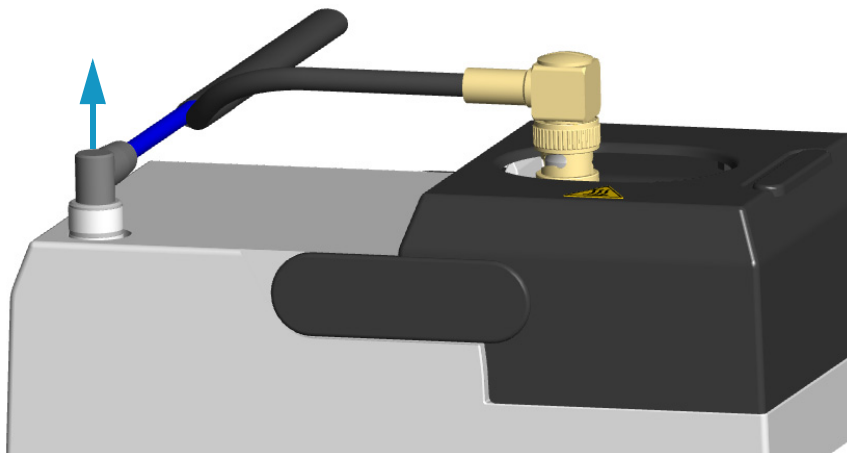


8. Open the module flap cover.
9. Remove the thermionic source assembly cable.
  - a. Twist the ring to disconnect the thermionic source assembly cable.
  - b. Push and twist the lock so that the button slides up in the groove, then pull the cable ends apart. See [Figure 127](#).

**Figure 127.** NPD Thermionic Source Cable Removal

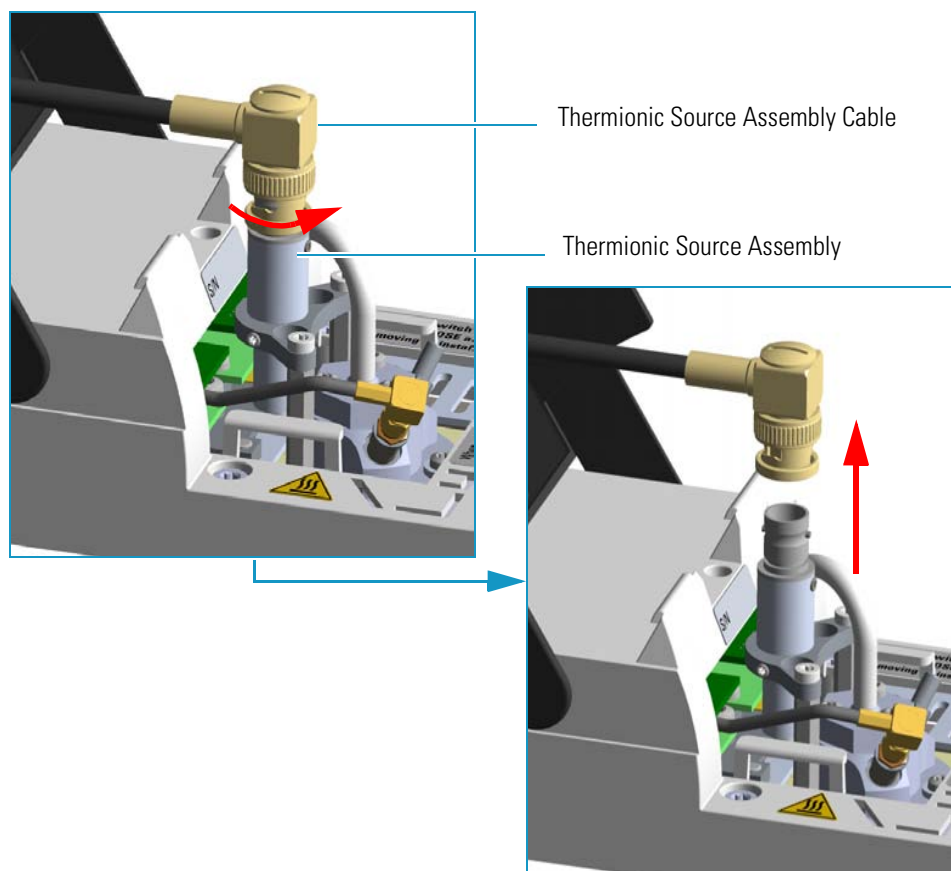


10. Unscrew and pull out the straight jack connector of the collecting electrode cable from the collecting electrode bulkhead jack. See [Figure 128](#).

**Figure 128.** Collecting Electrode Cable Removal

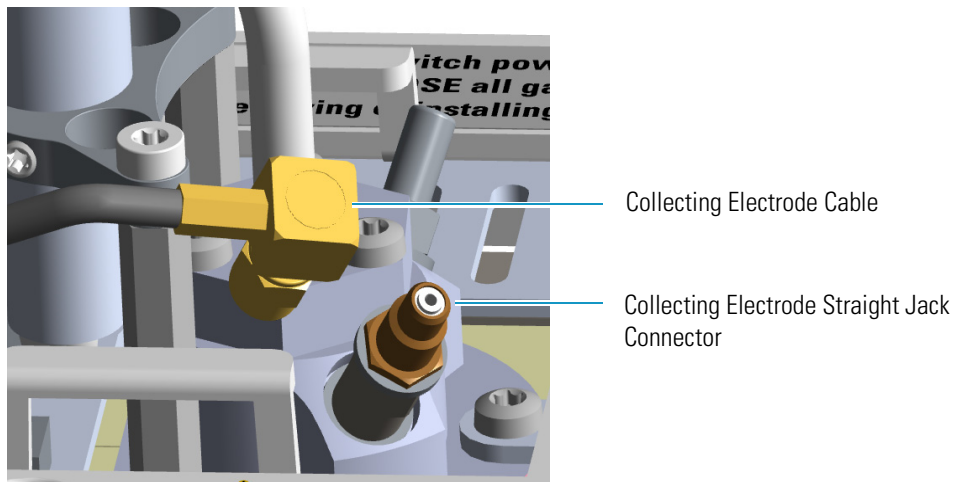


11. Remove the thermionic source.
  - a. Using a T6 Torxhead key, loosen the dowel which fix the thermionic source connector to the thermionic source support. See [Figure 129](#).

**Figure 129.** NPD Thermionic Source Removal (1)

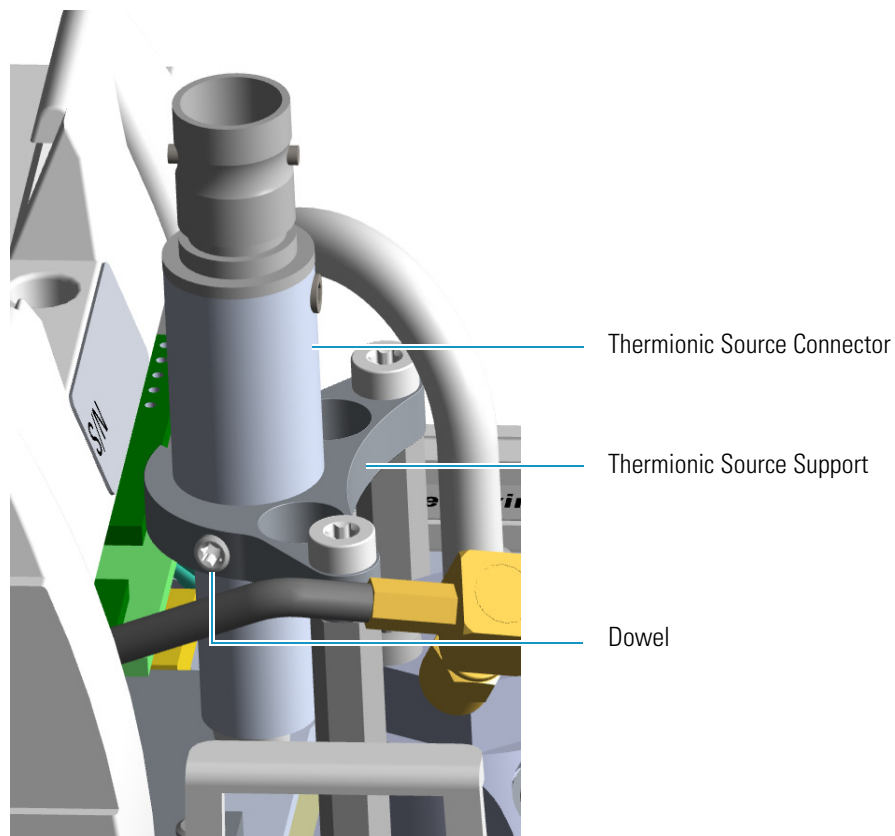


- b. Using the T10 Torxhead screwdriver, remove the three T10 Torx screws from the thermionic source assembly. See [Figure 130](#).

**Figure 130.** NPD Thermionic Source Removal (2)

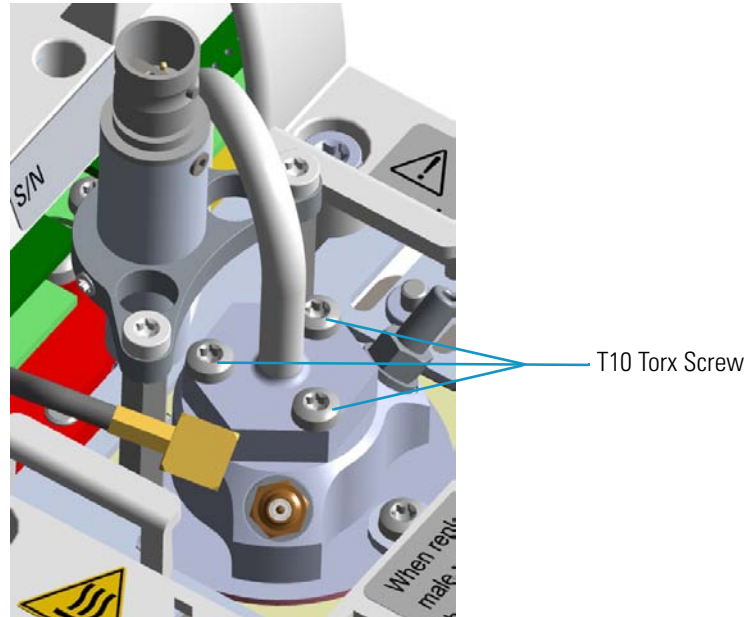

- c. Gently lift up the thermionic source connector guiding the flexible cable from the thermionic source support, then remove the thermionic source assembly from the NPD body. Avoid bumping the bead on the sides of the collector. See [Figure 131](#).

**Figure 131.** NPD Thermionic Source Removal (3)

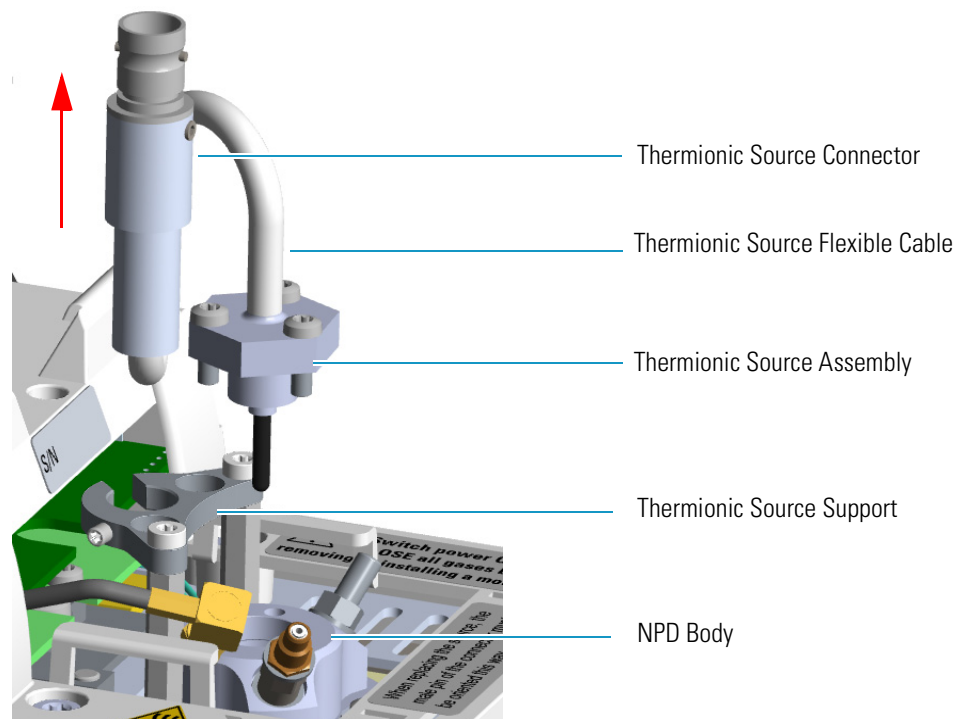


## 4 Performing Detectors Routine Maintenance

### Maintaining a Nitrogen Phosphorous Detector (NPD)


12. Replace the thermionic source.
  - a. Remove the protective cap covering the new thermionic source.
  - b. Mount and guide the new source assembly on the NPD body proceeding in the reverse order in which it was removed. Be careful not to bump the bead on the sides of the body and collecting electrode.
  - c. Rotate and align the thermionic source connector, then tighten the dowel to fix the connector to the support using the T6 Torxhead key. See [Figure 132](#).

**Figure 132.** Replace NPD Thermionic Source (1)

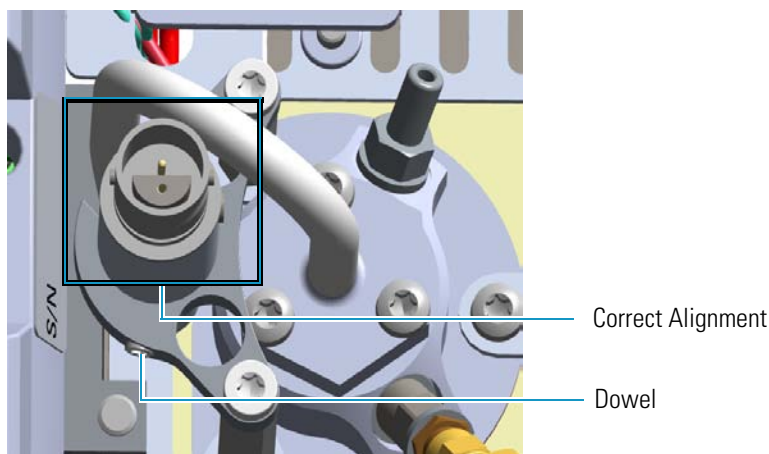

13. Reconnect and screw the straight jack connector of the collecting electrode cable to the collecting electrode bulkhead jack. See [Figure 133](#).


**Figure 133.** Replace NPD Thermionic Source (2)

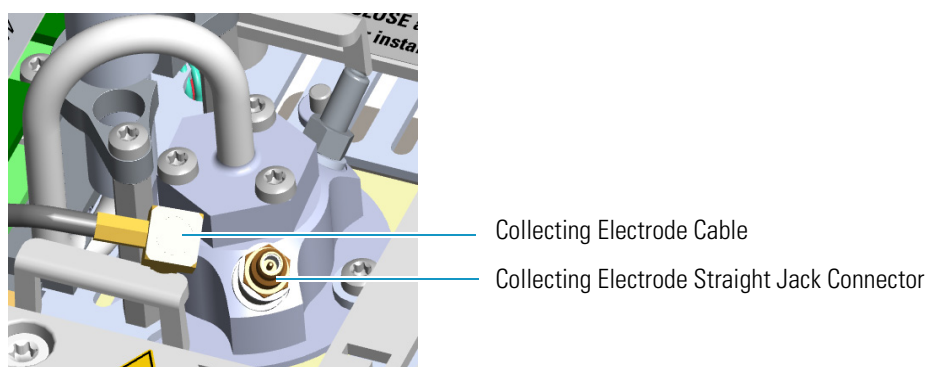

14. Reconnect the source assembly cable to the NPD thermionic source, and twist the ring to lock the connection. See [Figure 134](#).




**Figure 134.** Reconnect Source Assemble Cable

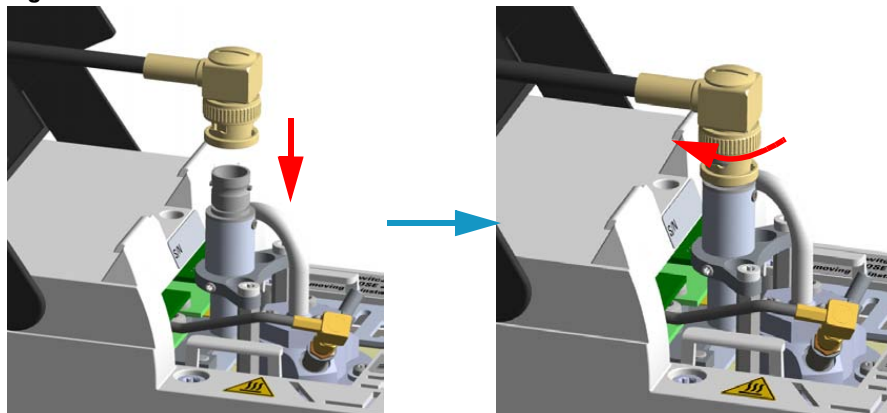

15. Close the module flap cover.
16. Plug in the signal cable into its contact on the detector module. See [Figure 135](#).

**Figure 135.** Reconnect Signal Cable

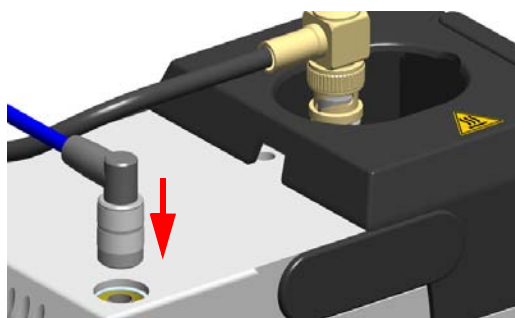

17. Plug the power cable to the AC input connector of the **NPD Thermionic Source Power Module** and to the wall outlet. The LED marked **On** blinks green.

If other external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.

18. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
19. Set the make-up gas on.
20. Switch on the thermionic source. Hydrogen and air are automatically opened. The LED marked **On** of the NPD Thermionic Source Power Module becomes solid green.
21. With all gases on, heat the detector to 150 °C and hold for about 15 minutes. Increase the temperature to 250 °C and hold for 15 minutes.

22. Increase the temperature to the operating value: 300 to 320 °C is recommended. Allow 15 minutes for equilibration.
23. Check the NPD leakage current. If > 2.0 pA, verify the installation.

## Cleaning or Replacing the NPD Collecting Electrode

### ❖ To clean or replace the NPD collecting electrode

1. Put the GC in standby condition.
2. Switch off the thermionic source. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.
4. Set the make-up gas off.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

5. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
6. Unplug the power cable from the AC input connector of the **NPD Thermionic Source Power Module** and from the wall outlet.

If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.

7. Unplug the signal cable from the detector module. See [Figure 136](#).


**Figure 136.** Unplug NPD Signal Cable

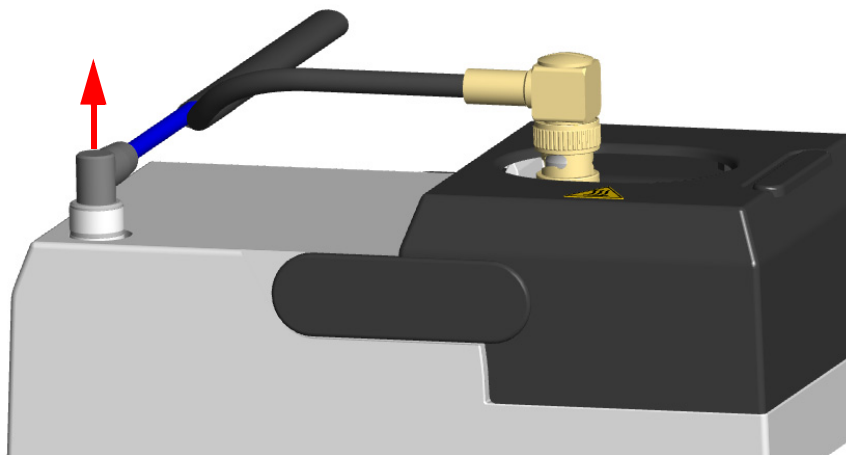

8. Open the module flap cover.


9. Remove the thermionic source assembly cable.
  - a. Twist the ring to disconnect the thermionic source assembly cable.
  - b. Push and twist the lock so that the button slides up in the groove, then pull the cable ends apart. See [Figure 137](#).

**Figure 137.** NPD Thermionic Source Cable Removal

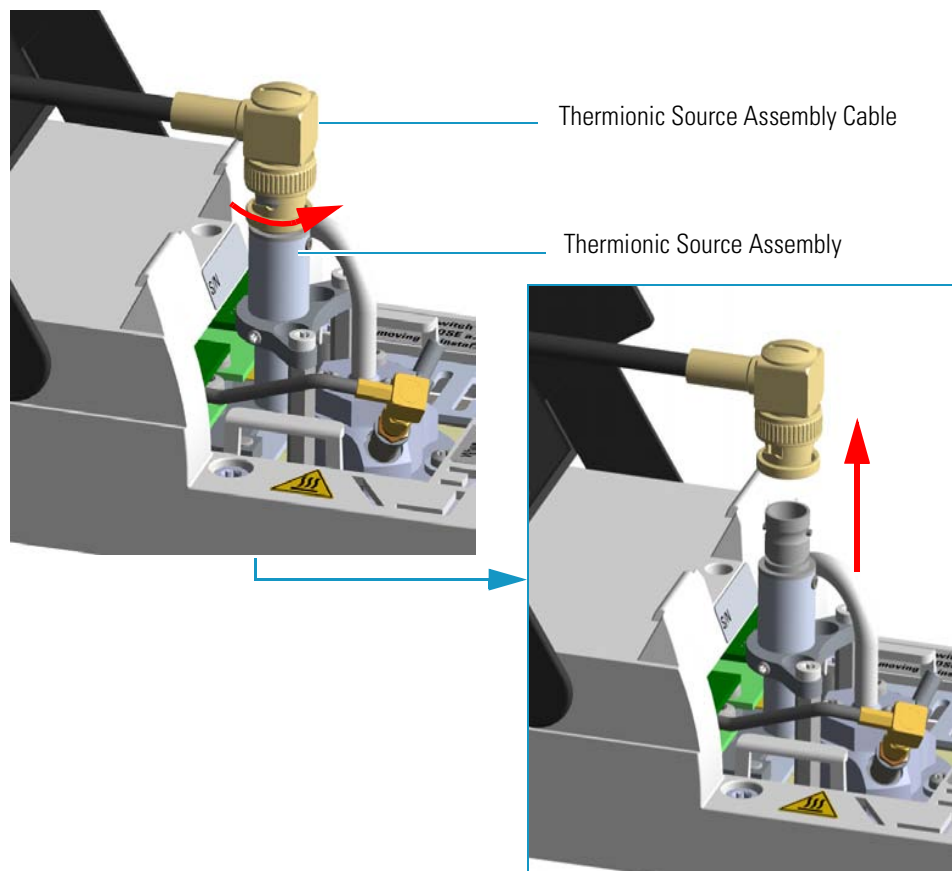

10. Unscrew and pull out the straight jack connector of the collecting electrode cable from the collecting electrode bulkhead jack. See [Figure 138](#).

**Figure 138.** Collecting Electrode Cable Removal

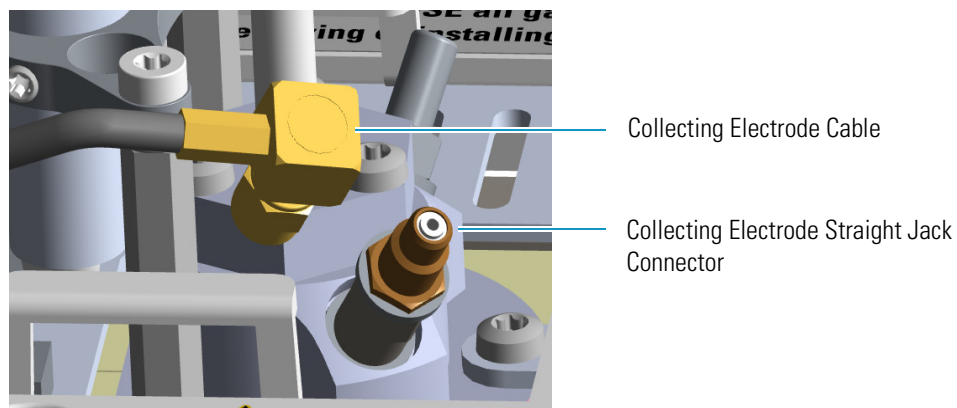

## 4 Performing Detectors Routine Maintenance

### Maintaining a Nitrogen Phosphorous Detector (NPD)

11. Remove the thermionic source.




**CAUTION** The thermionic source is delicate. Be careful not to break or crack the source. When performing maintenance on the NPD, avoid touching the source with your fingers, and prevent it from coming in contact with other surfaces.

Wear clean, lint-free gloves to prevent contamination of parts with dirt and skin oils.


- a. Using a T6 Torxhead key, loosen the dowel which fix the thermionic source connector to the thermionic source support. See [Figure 139](#).

**Figure 139.** NPD Thermionic Source Removal (1)

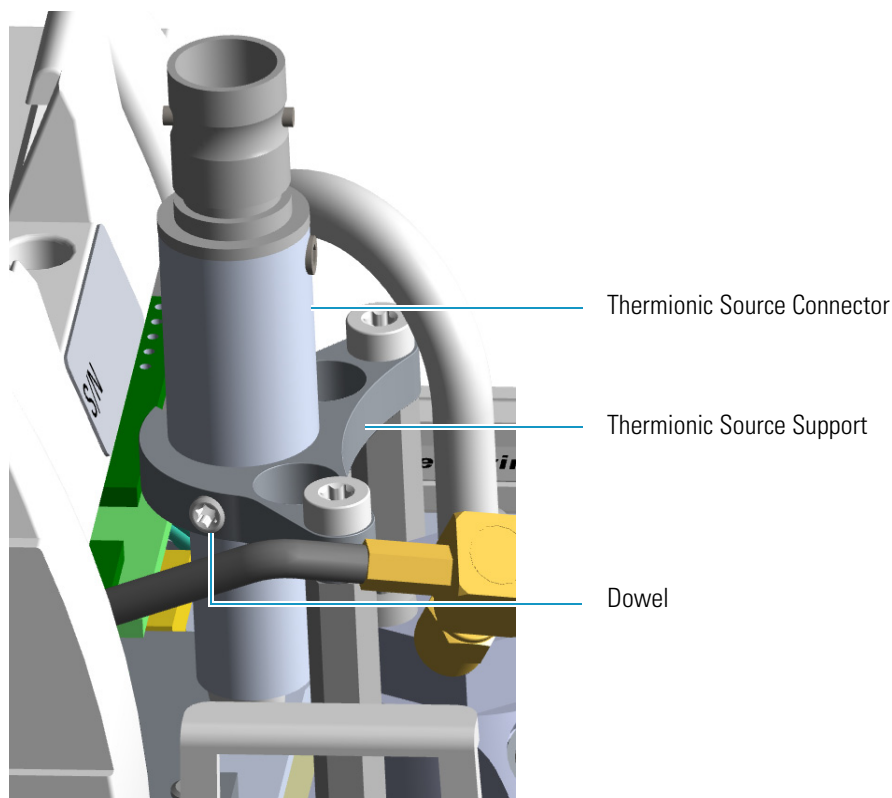


- b. Using the T10 Torxhead screwdriver, remove the three T10 Torx screws from the thermionic source assembly. See [Figure 140](#).

**Figure 140.** NPD Thermionic Source Removal (2)

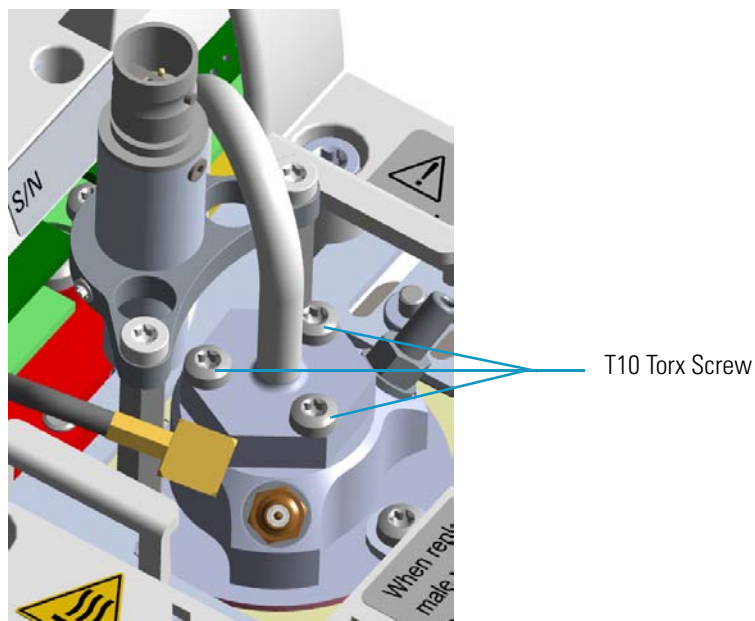

- c. Gently lift up the thermionic source connector guiding the flexible cable from the thermionic source support, then remove the thermionic source assembly from the NPD body. Avoid bumping the bead on the sides of the collector. See [Figure 141](#).

**Figure 141.** NPD Thermionic Source Removal (3)

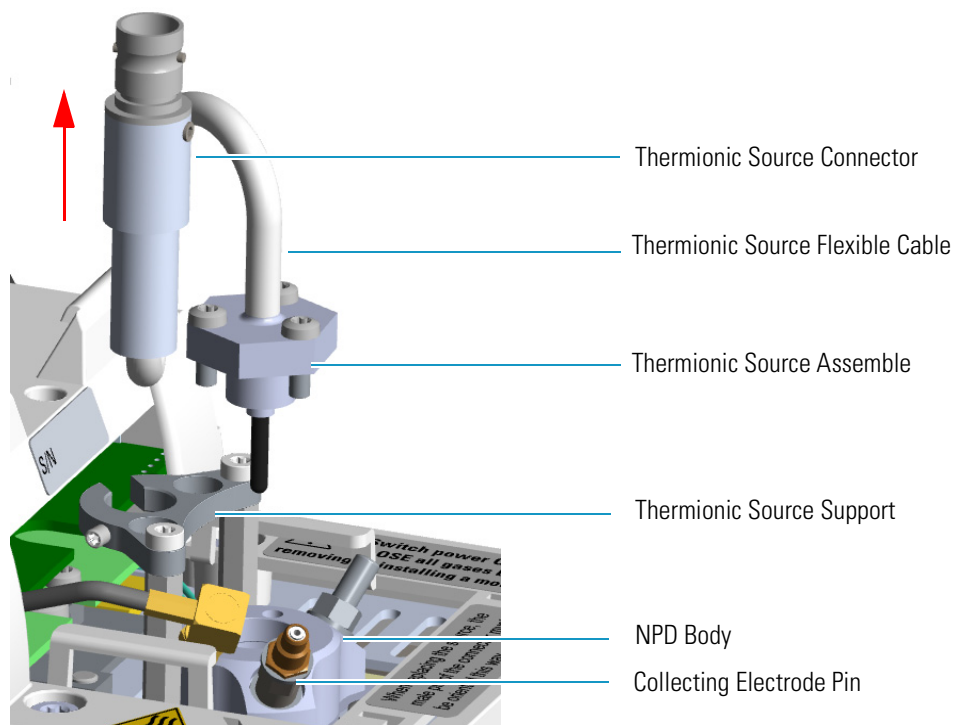


12. Remove and clean the collecting electrode.

## 4 Performing Detectors Routine Maintenance

### Maintaining a Nitrogen Phosphorous Detector (NPD)

- a. Using a 6-mm wrench, unscrew and remove the collecting electrode pin from its seat on the detector body.
- b. Using forceps or tweezers, extract the collecting electrode through the top of the detector body. See [Figure 142](#).

**Figure 142.** NPD Collecting Electrode Removal

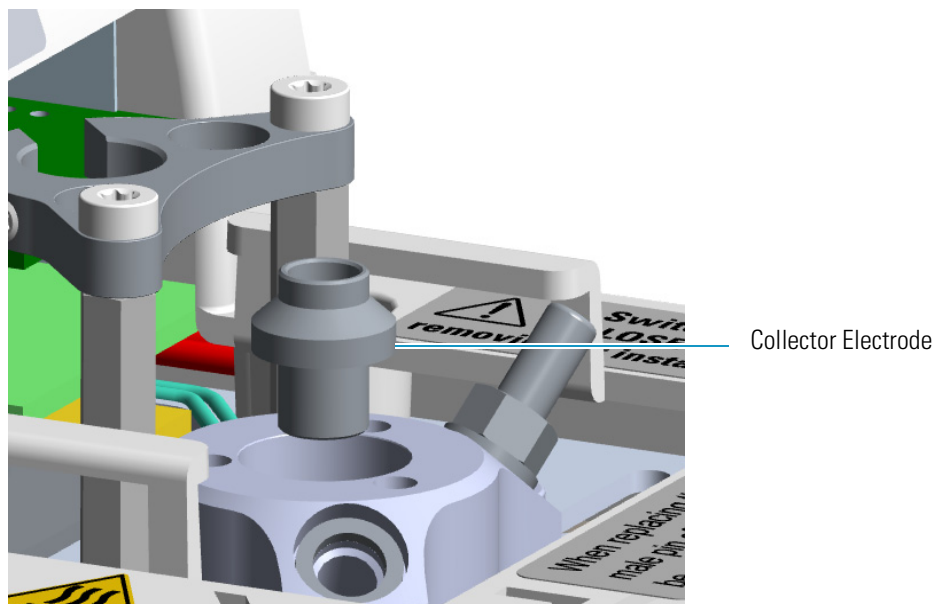


- c. Place the collecting electrode in the ultrasonic cleaner filled with liquid detergent.
- d. Sonicate the collecting electrode for about five minutes.
- e. Handling the collecting electrode with forceps or tweezers, rinse it using distilled water, then methanol.
- f. Place the collecting electrode on a paper towel and let it air dry.

**Note** If after cleaning the physical condition of the collecting electrode does not permit its reuse, replace the electrode with a new one.

13. Reinstall the collecting electrode.
  - a. Place the cleaned (or a new) collecting electrode in its housing.
  - b. Using a 6-mm wrench, screw the collecting electrode pin into its seat on the detector body.
  - c. Screw the straight plug crimp connector of the collecting electrode cable to the collecting electrode bulkhead jack.
14. Reinstall the thermionic source.
  - a. Remove the protective cap covering the new thermionic source.


- b. Mount and guide the new source assembly on the NPD body proceeding in the reverse order in which it was removed. Be careful not to bump the bead on the sides of the body and collecting electrode.
- c. Rotate and align the thermionic source connector, then tighten the dowel to fix the connector to the support using the T6 Torxhead key. See [Figure 143](#).

**Figure 143.** Replace NPD Thermionic Source (1)

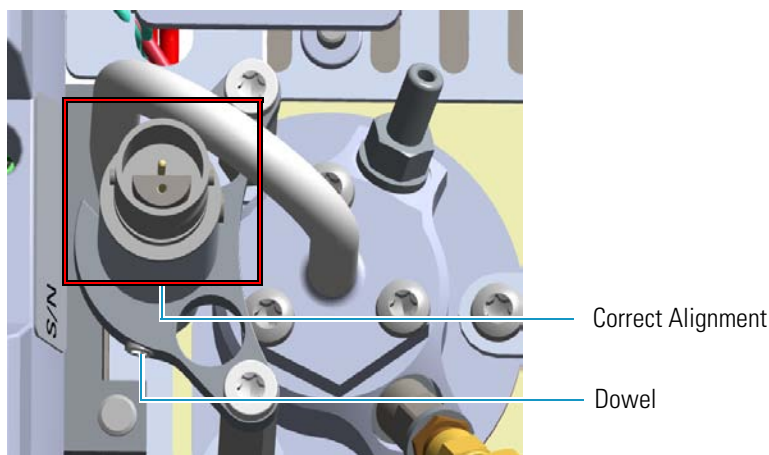

15. Reconnect and screw the straight jack connector of the collecting electrode cable to the collecting electrode bulkhead jack. See [Figure 144](#).


**Figure 144.** Replace NPD Thermionic Source (2)

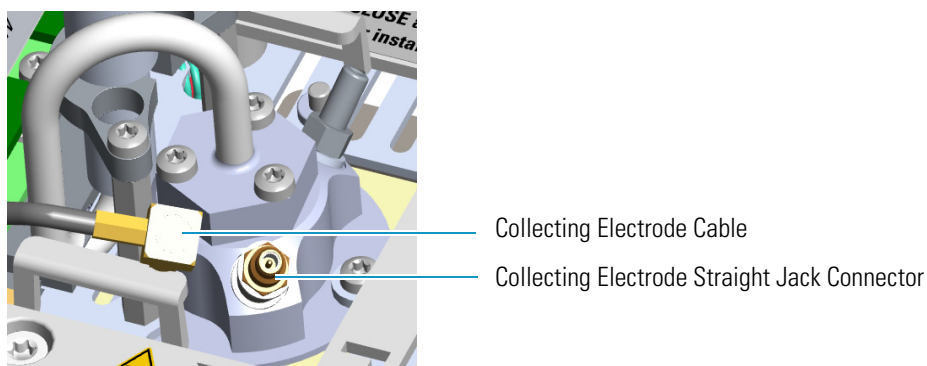

16. Reconnect the source assemble cable to the NPD cable and twist the ring to lock the connection. See [Figure 145](#).




**Figure 145.** Reconnect Source Assemble Cable

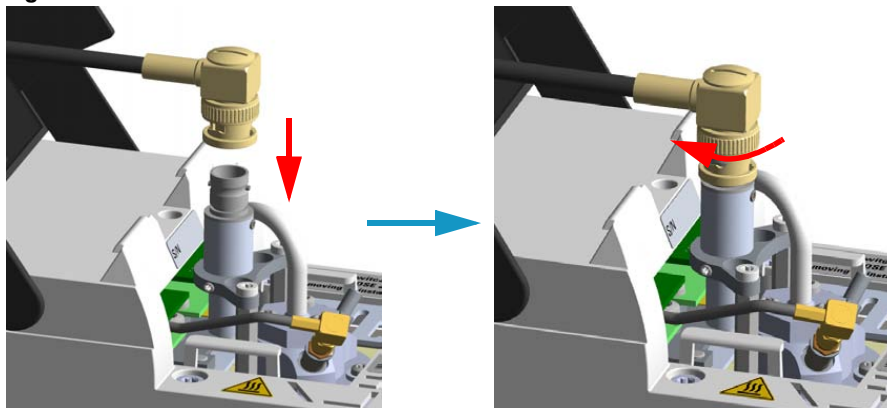

17. Close the module flap cover.
18. Plug in the signal cable into its contact on the detector module. See [Figure 146](#).

**Figure 146.** Reconnect Signal Cable

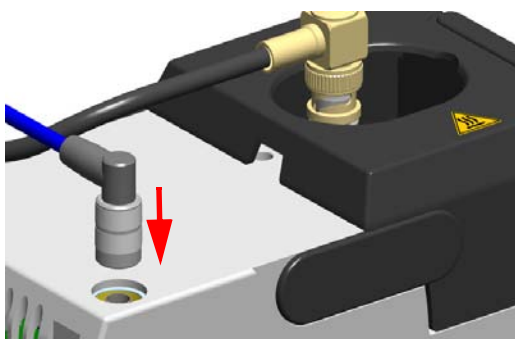


19. Plug the power cable to the AC input connector of the **NPD Thermionic Source Power Module** and to the wall outlet. The LED marked **On** blinks green.  
  
If other external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.
20. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
21. Turn the make-up gas on.
22. Heat the detector at the working temperature.
23. Switch on the thermionic source. Hydrogen and air are automatically opened.  
The LED marked **On** of the NPD Thermionic Source Power Module becomes solid green.



## Maintaining a Thermal Conductivity Detector (TCD)


This section provides instructions for maintaining a Thermal Conductivity Detector (TCD). The component of module and detector are shown in [Figure 147](#).

**Figure 147.** TCD Module Components

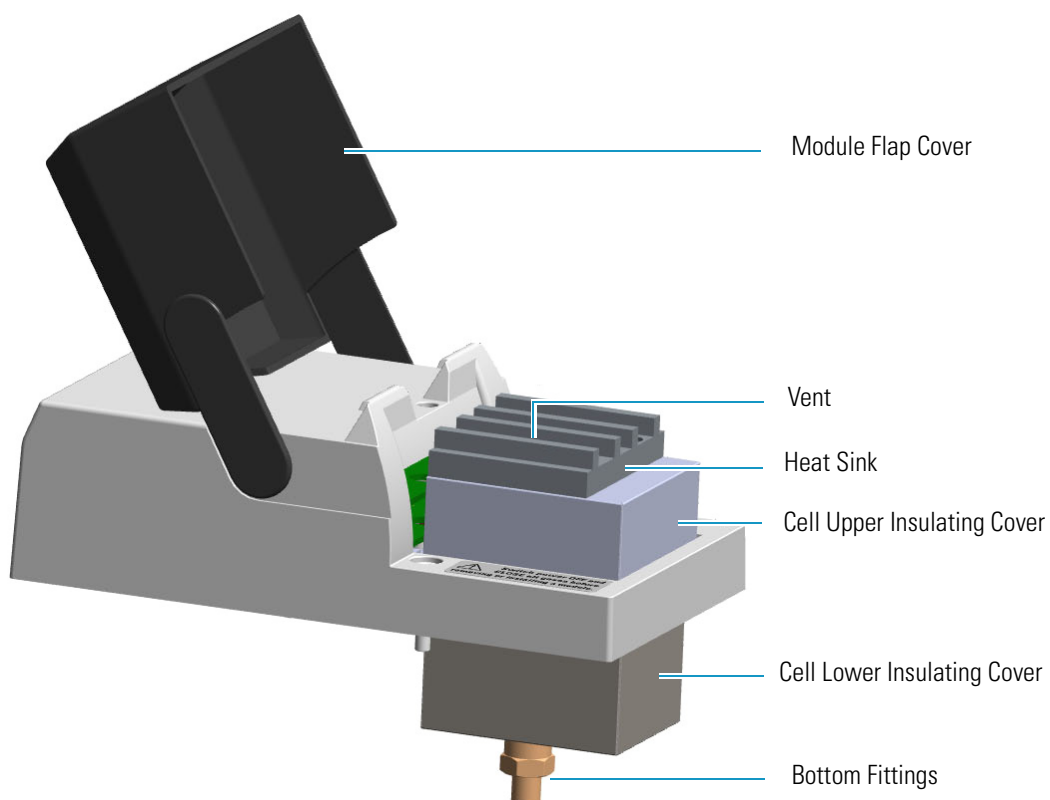

The same routine maintenance is applied to the **In-Series Connection TCD Module**. See [Figure 148](#).

**Figure 148.** In-Series Connection TCD Module

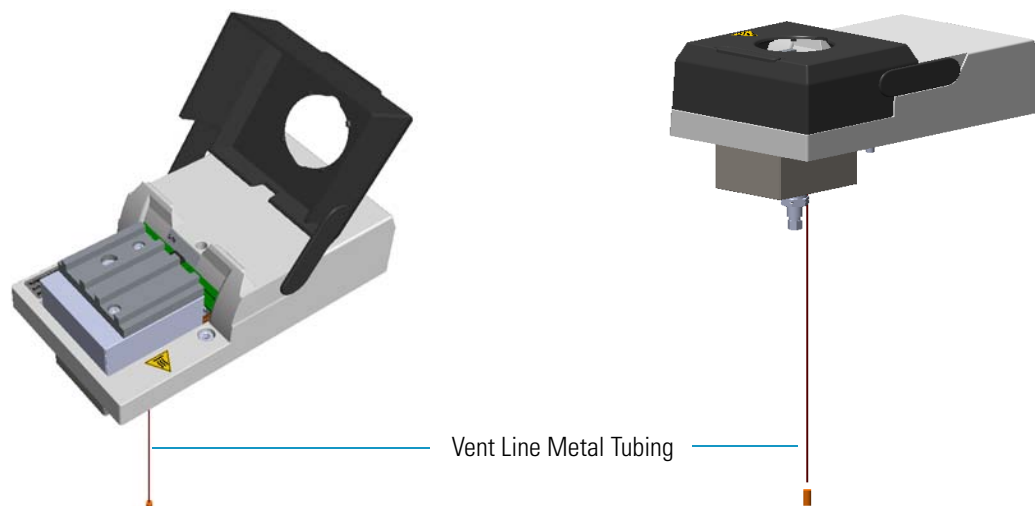

## 4 Performing Detectors Routine Maintenance

### Maintaining a Thermal Conductivity Detector (TCD)

The TCD detector does not usually need maintenance. Nevertheless, if you follow a few simple precautions, you will avoid problems and prolong the detector's lifetime. Pay a special attention to avoid contamination or damage to the filaments.

The detector must be returned to the factory when filament replacement is required.



**CAUTION** Do not turn the filament on until the carrier gas is flowing through the detector.

Do not shut off or disconnect the carrier gas when the detector is hot, even if the unit is turned off.

Follow these simple rules:

- Avoid switching on the filaments when unnecessary. Doing so considerably reduces the filament's lifetime.
- Avoid injecting samples that contain halogenated or acid compounds at high concentrations.
- Ensure that oxygen (air) cannot enter into the filament's cells. Oxidation will irreversibly damage the filaments. Install traps for moisture and oxygen on the gas lines to reduce the risk of contamination.
- Turn off the filament before disconnecting the column from the detector. When the column is disconnected, air will enter into the cell and the filament, if powered on, will be damaged.
- Turn on the filament only if the column is connected. It is a good practice to let the gas flow through the cells for 5-10 minutes before powering the filaments.

See [Bake-out Procedure](#), [Measuring the Carrier Gas Flow Rate](#), and [Shutting Down the TCD](#) for additional information.

## Bake-out Procedure

Under normal conditions, the TCD requires no routine maintenance. However, if the detector is exposed to chemicals which may condense or polymerize within the detector and adversely affect performance, an attempt of filament cleaning could be baked out at high temperatures (up to 300°C). Inert carrier gas flow should be maintained during the reconditioning procedure.

Filament temperature should also be increased to approximate the setpoint of the bake-out temperature.

If a 24-hour bake-out is not sufficient to remove the contamination, the unit must be returned to the factory for disassembly and cleaning.


## Measuring the Carrier Gas Flow Rate

While the measuring of the carrier gas flow rate through the **standard TCD detector** is performed simply connecting the flow-meter to the vent outlet, the **TCD In-series detector** requires the use of the **column- flowmeter connector**.

### ❖ To measure the TCD In-Series carrier gas flow rate

Carefully push the vent line metal tubing end into the flow meter section of the column-flowmeter connector. See [Figure 149](#).

**Figure 149.** Column Flowmeter Connector

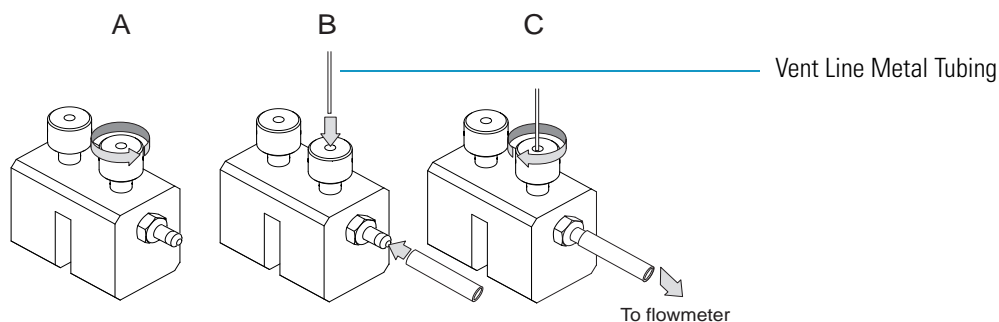

Connect the flowmeter to the dedicated fitting on the column-flowmeter connector, then measure the carrier gas flow rate.

## Shutting Down the TCD


To shut down the detector at the end of the analytical cycle:

1. Turn off the filaments.
2. Reduce the carrier gas flow to 50% of the normal operating flow to conserve gas supplies.


## Maintaining an Electron Capture Detector (ECD)

This section provides instructions for maintaining an Electron Capture Detector (ECD). The component of module and detector are shown in [Figure 150](#) and [Figure 151](#).

**Figure 150.** ECD Module Components

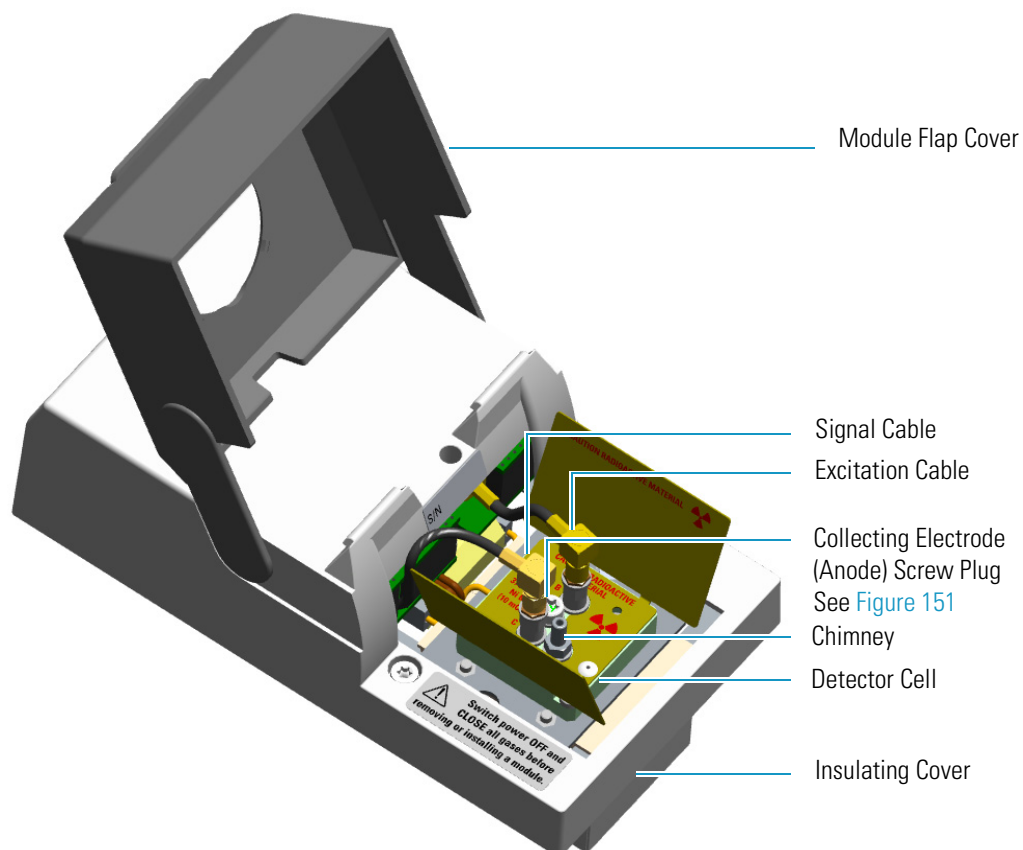

**Figure 151.** Collecting Electrode (Anode) Screw Plug

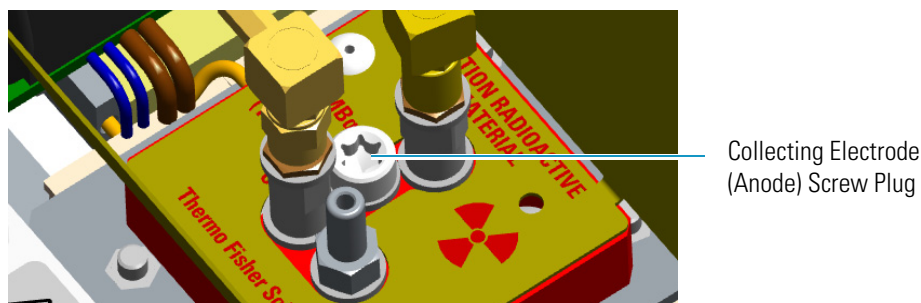

The ECD detector periodic maintenance includes:

- **Cleaning the collecting electrode (Anode)**

We suggest cleaning the collecting electrode annually or when a high noise baseline is found producing detector contamination. Replace the collecting electrode if necessary. See [“Cleaning or Replacing the Collecting Electrode \(Anode\)”](#) on [page 169](#).

Before maintaining the detector, read the following precautions and notes:


**CAUTION** The Electron Capture Detector contains a  $^{63}\text{Ni}$  beta-emitting radioactive source at 370 MBq (10 mCi). The detector should never be opened or handled by the user. Any maintenance or service operations involving even partial disassembling of the detector must be performed **ONLY** by qualified personnel at a laboratory expressly authorized by Thermo Fisher Scientific and specifically licensed to handle radioactive material.

**IMPORTANT** For customers within the jurisdiction of the United States Nuclear Regulatory Commission (US NRC), you may find a listing of the agreement states and the current contact information for the regulators covering both Generally and Specifically licensed devices like an ECD at <http://nrc-stp.ornl.gov/rulemaking.html>. This information is maintained by the US NRC.


**IMPORTANT** For customers within the jurisdiction of the Canadian Nuclear Safety Commission, you may find current contact information for the regulators at <http://www.nuclearsafety.gc.ca>

**IMPORTANT** The recommended working life of the detector is 15 years, after which the user should arrange for the radioactive source to be inspected and assessed by a qualified authority to extend its working life, or dispose of the product through a suitable disposal route.

**IMPORTANT** According to International, USA, and Canada licenses, the appropriate Radioactivity Data labels are applied on the ECD detector module.

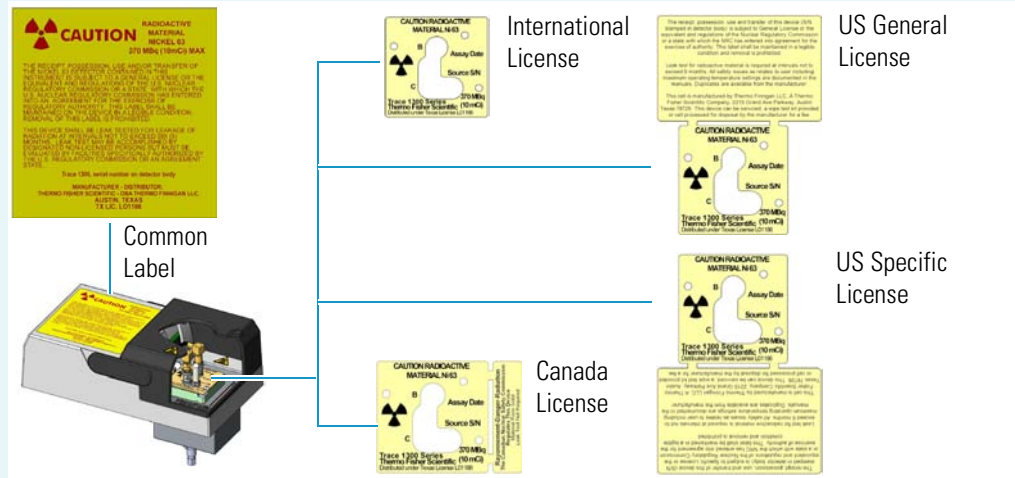

Exclusively for ECD detector module used in **China**, the **Caution Radioactive Material Ni63** labels are replaced by the followings:

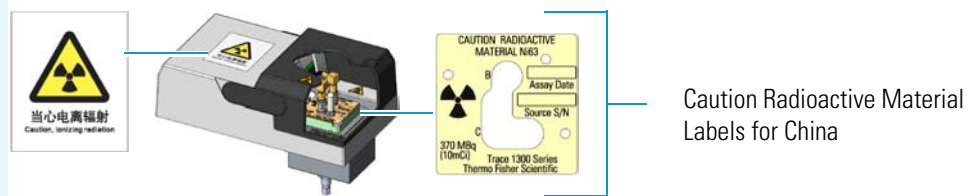

## 4 Performing Detectors Routine Maintenance

### Maintaining an Electron Capture Detector (ECD)



**WARNING** All the operations must be carried out at low temperature to avoid burns. Therefore, before conducting maintenance, cool the detector to room temperature.

#### Materials needed to maintain an ECD detector

Ultrasonic cleaner

GC-grade hexane

GC-grade toluene

Fine emery paper

Little flathead screwdriver

T20 Torxhead screwdriver

Forceps or tweezers

Collecting electrode (Anode) if necessary

## Detector Chemical Contamination

The ECD, if properly used, has a good resistance against chemical contamination. However, some critical operating conditions may cause, over time, contamination of the collecting electrode (anode). This contamination is highlighted by an excessive increase of the base frequency, and a baseline drift when the reference current or the pulse amplitude is changed. The collecting electrode can be easily removed and cleaned without disturbing the radioactive source.

Detector contamination could be indicated by the following effects in the chromatogram:

- Reduced signal to noise ratio
- High-noise baseline (high frequency)
- Baseline drift with changing pulse voltage
- Negative dips after peaks

If chemical contamination of the whole cell is suspected proceed as follows:

1. Heat the ECD at the maximum operating temperature with carrier and make-up gases flowing through the detector (thermal cleaning).
2. Follows the decontamination process by monitoring the base frequency. Initially the frequency value tends to increase to very high values, and then it decreases to acceptable ones.

If irreversible contamination of the cell is suspected please contact your local Thermo Fisher Scientific Technical Service office.

## Wipe Test

Before leaving the factory, each ECD is leak tested for surface radio contamination by means of the **Wipe Test** method. Each detector is provided with a **Wipe Test Certificate** reporting the results of the values found and the procedure followed.



**IMPORTANT** The users of the ECD detector in the United States are required to perform a **Wipe Test** on their ECD at intervals not exceeding 6 months, following the reported procedure. For other countries, please refer to the appropriate agency for equivalent requirements.

## Cleaning or Replacing the Collecting Electrode (Anode)

### ❖ To clean or replace the collecting electrode

1. Put the GC in standby condition.
2. Cool the detector to room temperature.
3. Turn the make-up gas off.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
6. Open the module flap cover.
7. Disconnects the signal and excitation cables from their contacts on the cell top cover.  
See [Figure 152](#).

## 4 Performing Detectors Routine Maintenance

Maintaining an Electron Capture Detector (ECD)


**Figure 152.** ECD Cables Disconnection

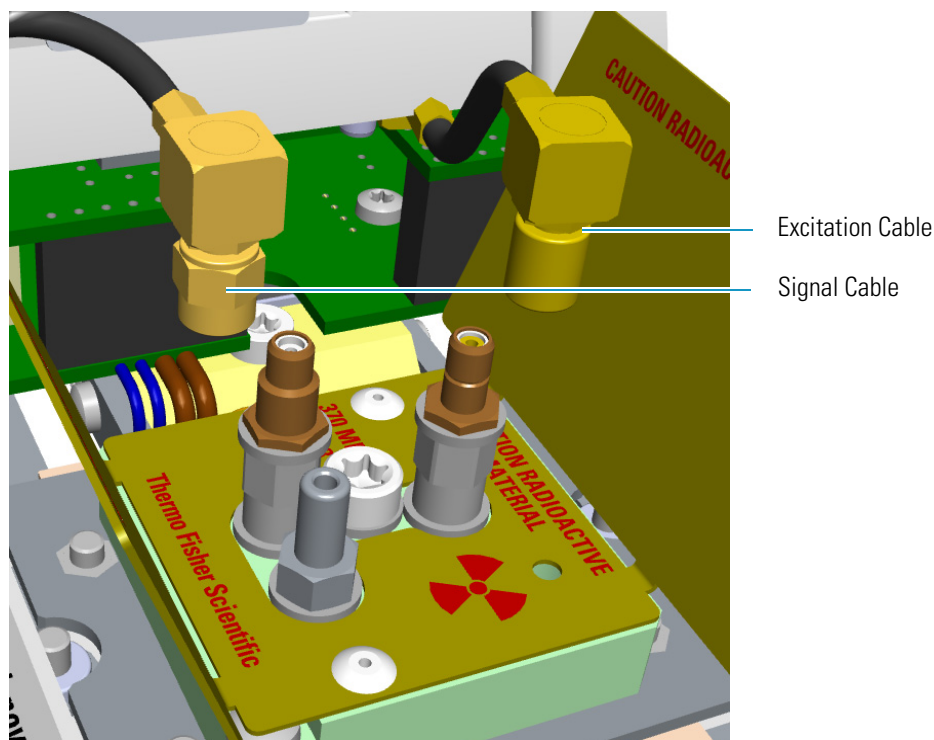

- a. Unscrew and pull out the straight plug crimp connector of the signal cable from the collecting electrode (anode) bulkhead jack.
  - b. Pull out the quick coupling straight jack connector of the excitation cable from the polarizing electrode bulkhead jack.
  - c. Carefully move the cables in order to have free space for handling the detector.
8. Remove and clean the collecting electrode (anode).
- a. Unscrew the screw plug with its seal hiding the collecting electrode (anode) using the T20 Torxhead screwdriver.
  - b. Using a little flathead screwdriver, unscrew and remove the collecting electrode (anode). See [Figure 153](#).



**Figure 153.** Collecting Electrode (Anode) Removal

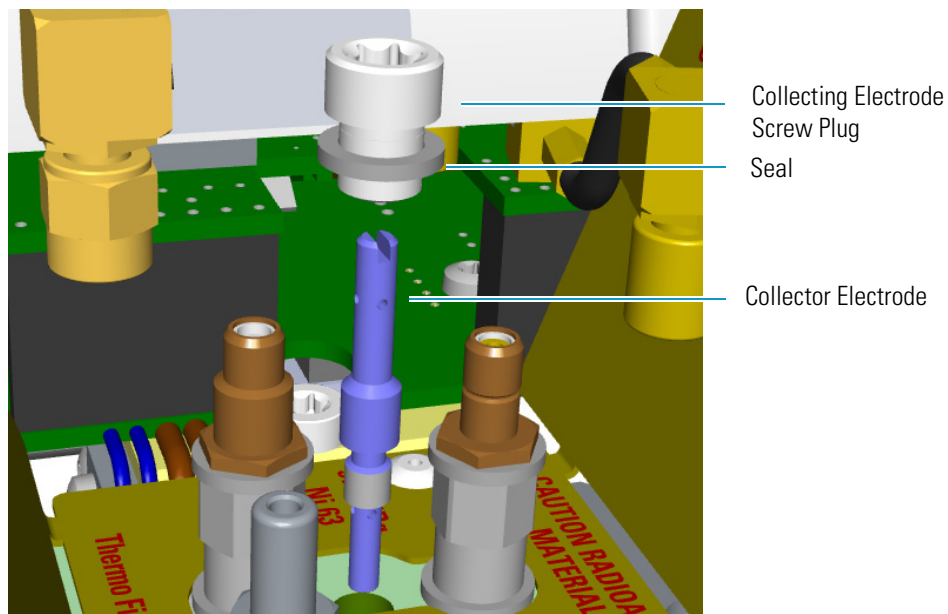


- c. Place the collecting electrode in the ultrasonic cleaner filled with liquid detergent, and clean it for about five minutes.
  - d. Handling the collecting electrode (anode) with forceps or tweezers, rinse it using distilled water, then methanol.
  - e. Place the collecting electrode (anode) on a paper towel and let it air dry.
- Note** If after cleaning the physical condition of the collecting electrode (anode) does not permit its reuse, replace the electrode with a new one.
9. Reinstall the collecting electrode (anode).
    - a. Screw the cleaned, or a new collecting electrode (anode) into its housing. Screw the hiding screw plug and its seal.
  10. Reconnect signal and excitation cables.
  11. Close the module flap cover.
  12. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
  13. Power on the GC.
    - a. Plug the power cable to the AC Input connector on the back of the GC, and to the wall outlet.
    - b. Flip up the power switch (breaker) to the position I.
  14. Set the make-up gas on.
  15. Set the normal detector working conditions.

## Maintaining a Flame Photometric Detector (FPD)


This section provides instructions for maintaining a Flame Photometric Detector (FPD).  
See [Figure 154](#).

**Figure 154.** FPD Control Module

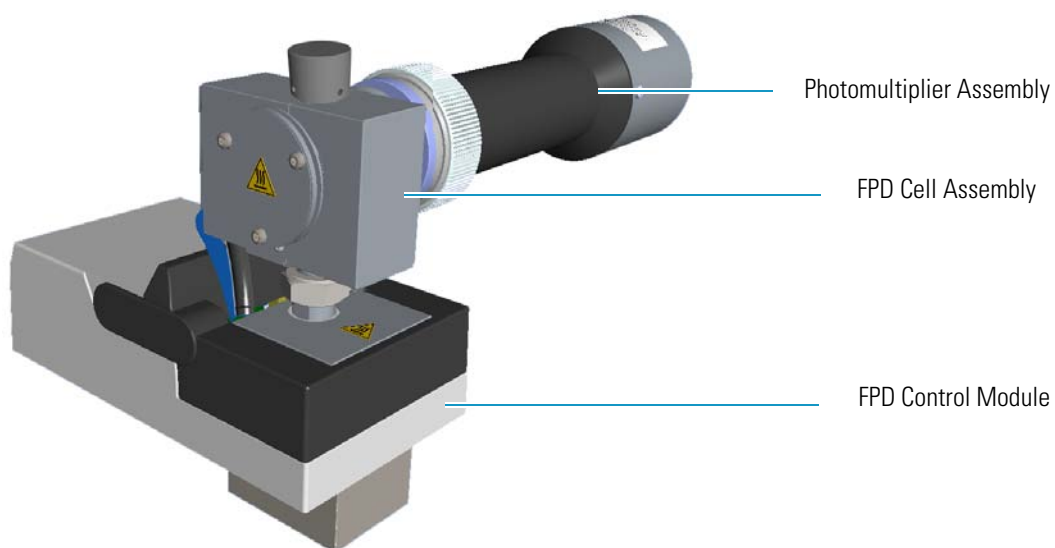


The control module and detector cell assembly components are shown in [Figure 155](#), [Figure 156](#), and [Figure 157](#).

**Figure 155.** FPD Control Module Components

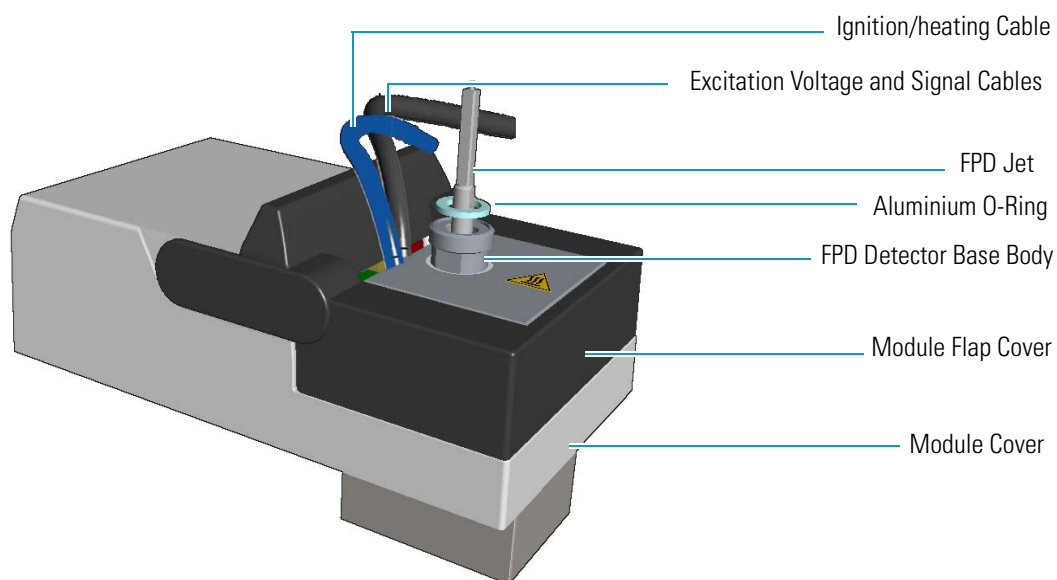

**Note** The signal, excitation voltage, and ignition/heating cables will be not shown in the next illustrations for graphic convenience.


**Figure 156.** FPD Detector Assembly Components

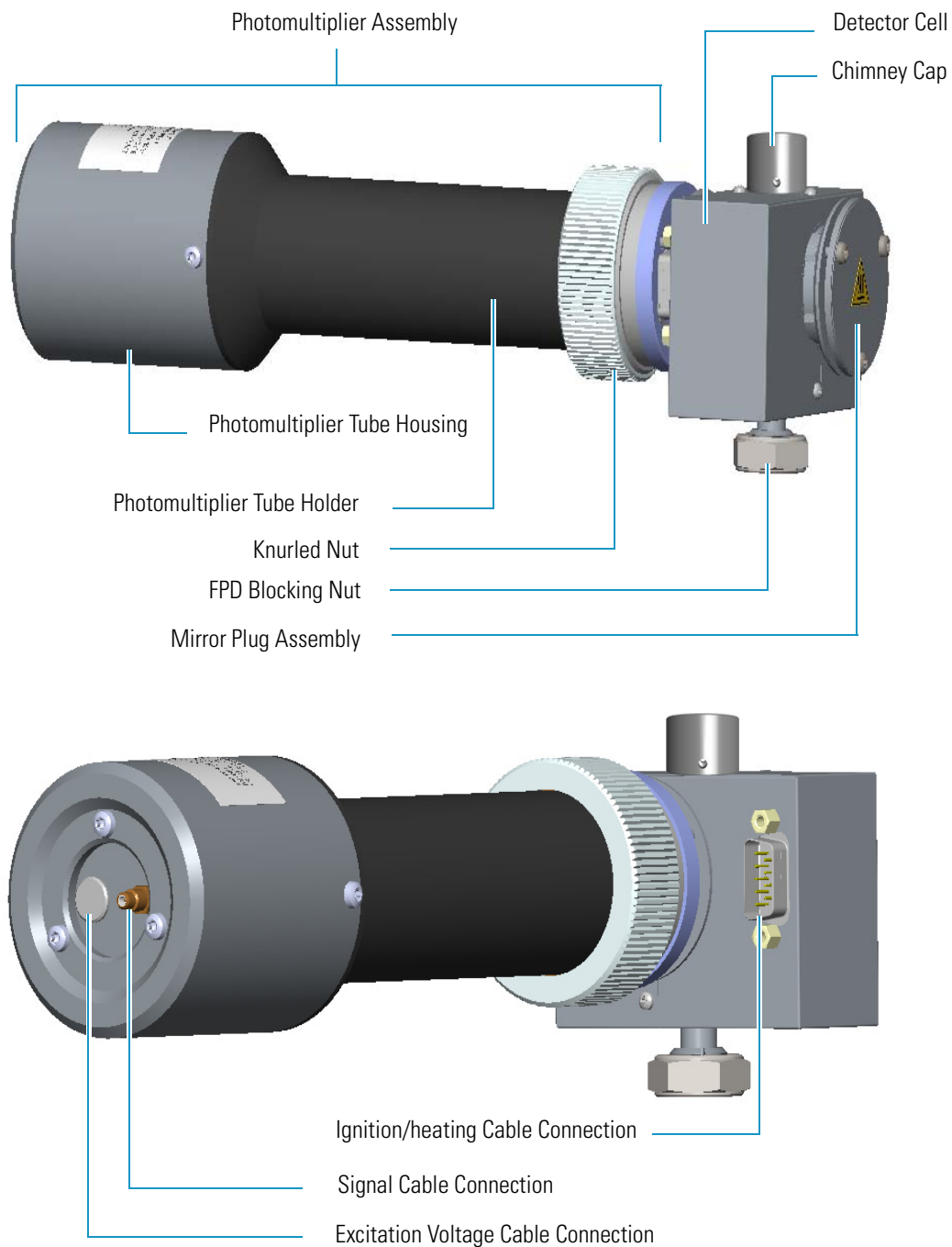

## 4 Performing Detectors Routine Maintenance

Maintaining a Flame Photometric Detector (FPD)

**Figure 157.** FPD Cell Assembly Components

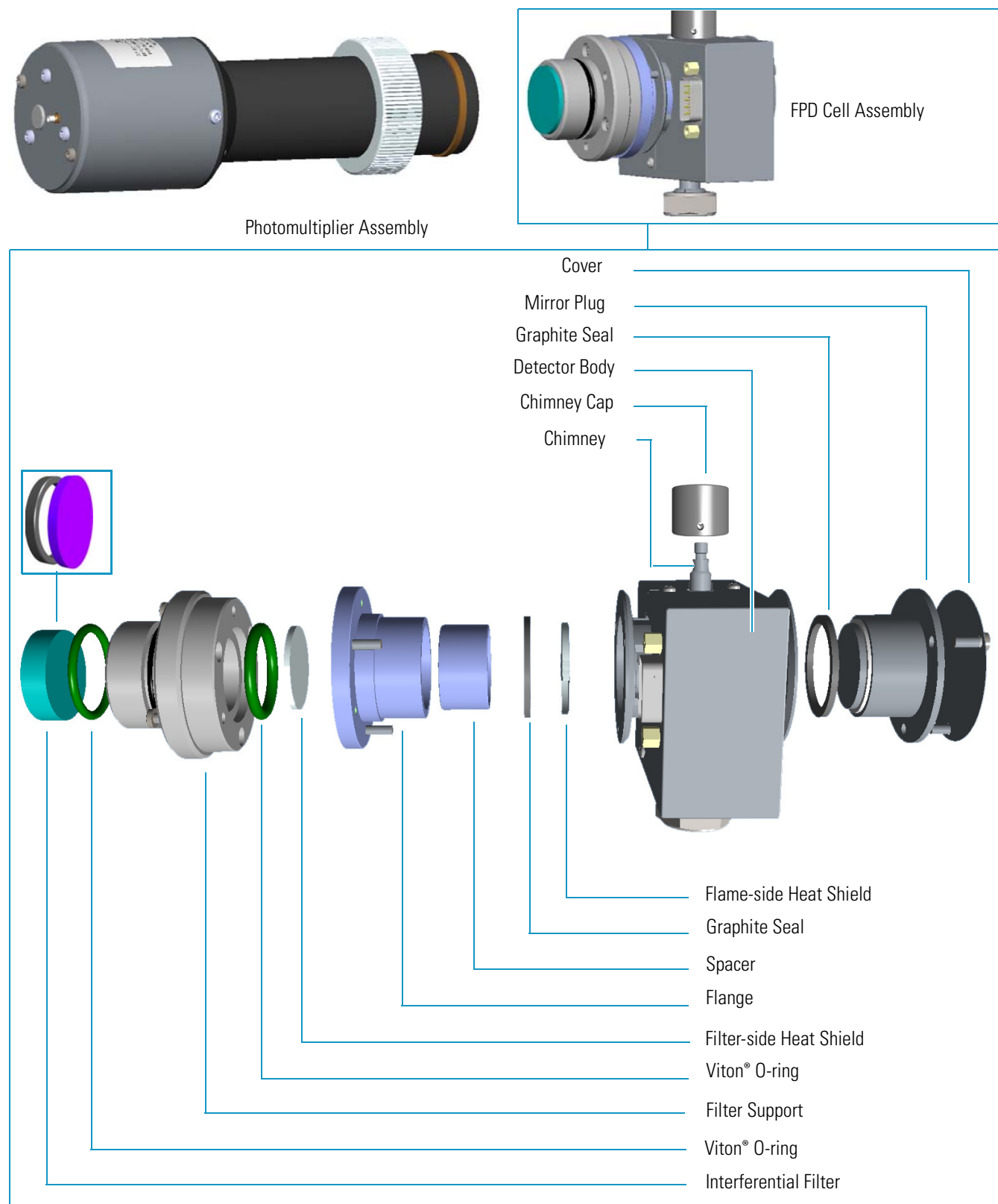

To ensure optimum performance of the FPD, you must keep it clean and free of dust and deposits. Symptoms such as reduced sensitivity and increased noise indicate that detector cleaning could be necessary.

The FPD detector periodic maintenance includes:

- **Installing the FPD detector**

See “[Installing the FPD Detector](#)” on [page 176](#).

- **Removing the FPD detector**

See “[Removing the FPD Detector](#)” on [page 178](#).

- **Cleaning or replacing the FPD jet**

You should clean the jet semiannually when you analyze particularly dirty compounds that the flame does not burn properly.

See “[Cleaning or Replacing the FPD Jet](#)” on [page 179](#).

- **Cleaning or replacing the FPD interferential filter**

See “[Cleaning or Replacing the FPD Interferential Filter](#)” on [page 183](#).

- **Replacing the FPD ignition glow-plug**

This operation is NOT part on the ordinary maintenance. The flame ignition element must be replaced only when defective.

See “[Replacing the FPD Ignition Glow-plug](#)” on [page 187](#).

Before maintaining the detector, read the following precautions:



**WARNING** Carry out all the operations at low temperature to avoid burns. Therefore, before beginning maintenance, cool the detector to room temperature.



**CAUTION** When handling organic solvents you must take precautions to avoid health hazards.

#### **Materials needed to maintain a FPD detector**

---

FPD fixing tool

---

Gloves

---

Ultrasonic cleaner

---

Liquid detergent

---

GC-grade methanol

---

Distilled water

---

## 4 Performing Detectors Routine Maintenance

### Maintaining a Flame Photometric Detector (FPD)

#### Materials needed to maintain a FPD detector

Paper towels

Forceps or tweezers

FPD Jet (if necessary)

Interferential filter for sulphur PN 28107000 or PN 19050785, (if necessary)

Interferential filter for phosphorous PN 28107100 or PN 19050785, (if necessary)

Interferential filter for tin PN 28107001 (optional), (if necessary)

FPD Ignition Glow-plug

T6 Torxhead key

T10 Torxhead screwdriver

5 mm wrench

5/16-in tube wrench

**Note** For maintaining or replacing any other component not listed in this section, see Chapter 7, “Detectors Advanced Maintenance.”

- “Measuring the FPD Gas Flows” on page 293
- “Cleaning or Replacing the FPD Mirror Metal Plug” on page 295
- “Cleaning or Replacing the FPD Filter-side Heat Shields” on page 300
- “Cleaning or Replacing the FPD Flame-side Heat Shields” on page 308
- “Replacing the FPD Photomultiplier Tube” on page 317


## Installing the FPD Detector

Before proceeding the installation of the FPD detector assembly, make sure that the FPD control module is properly installed into the back detector housing.

#### ❖ To install the FPD detector assembly on the FPD detector base body


1. Place the jet into the detector base body housing and tighten it. Ensure the jet is perfectly vertically aligned to avoid damage. [Figure 158](#).
2. Place the FPD on the detector base body, ensuring that the aluminium ring has been inserted in the correct position. Tighten the fixing nut by using the FPD fixing tool. See [Figure 158](#).

**Figure 158.** FPD Installation

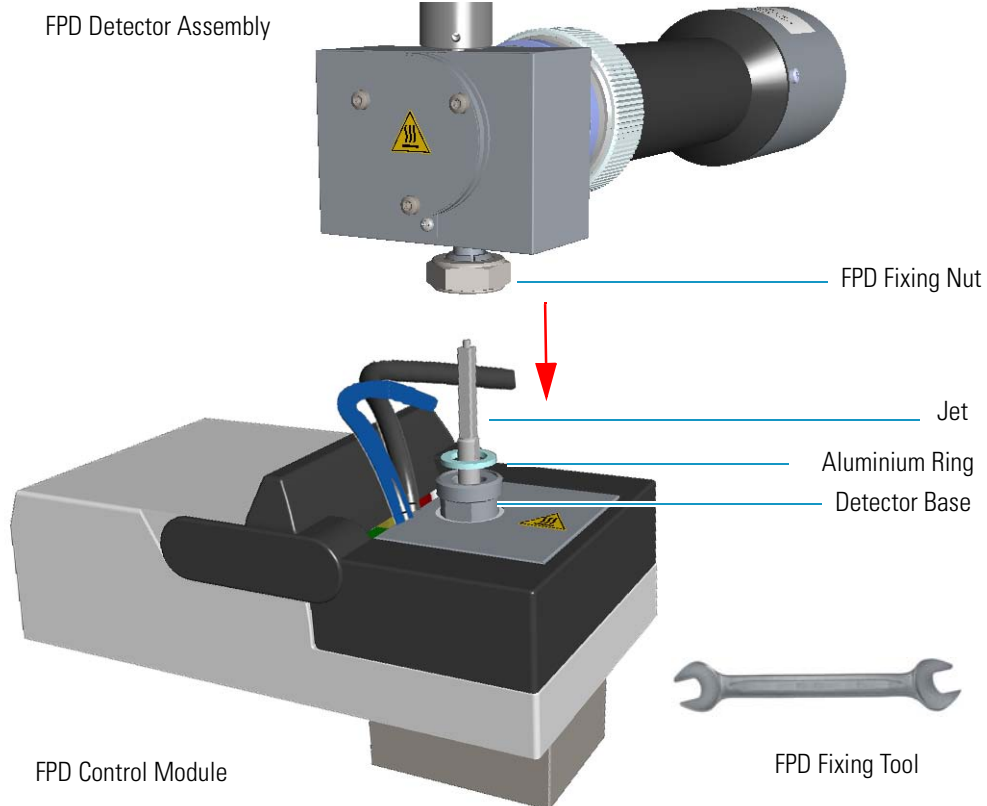

3. Carefully, connect the signal, excitation voltage, and ignition/heating cables, coming from the FPD control module, to the FPD detector. See [Figure 159](#).

**Figure 159.** FPD Cables Connection

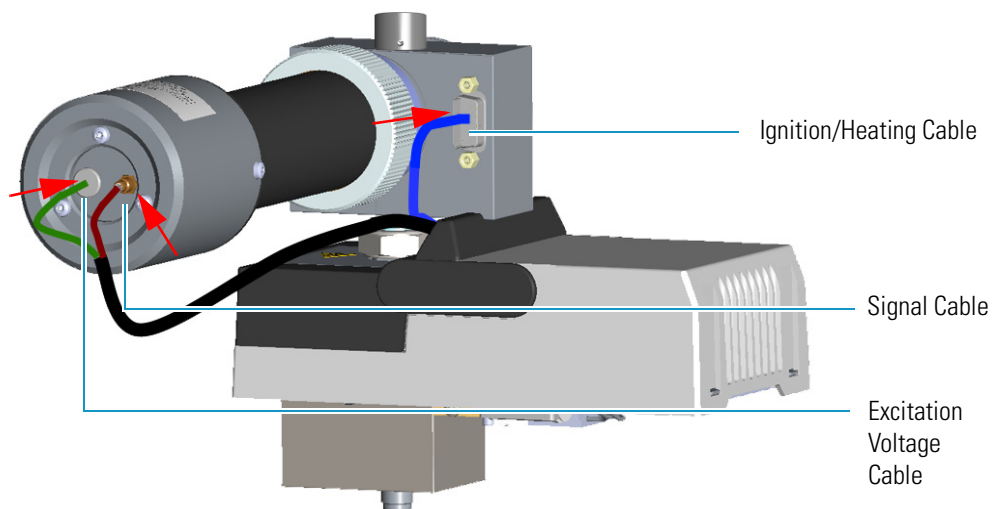

## Removing the FPD Detector


❖ **To remove the FPD detector assembly from the FPD detector base body**

1. Put the GC in standby condition.
2. Switch off the flame. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.
6. Remove the FPD detector from the base.
  - a. Disconnect the signal, excitation voltage, and ignition/heating cables from the detector. See [Figure 160](#)


**Figure 160.** FPD Cable Disconnection

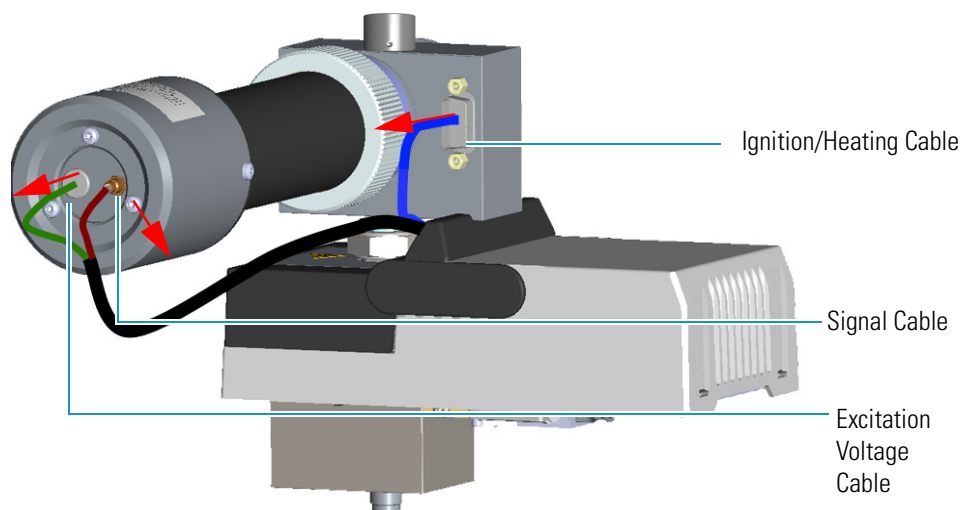

- b. Using the tool provided with the system, loosen the fixing nut on the base of the detector and remove it. See [Figure 161](#).



**Figure 161.** FPD Detector Removal

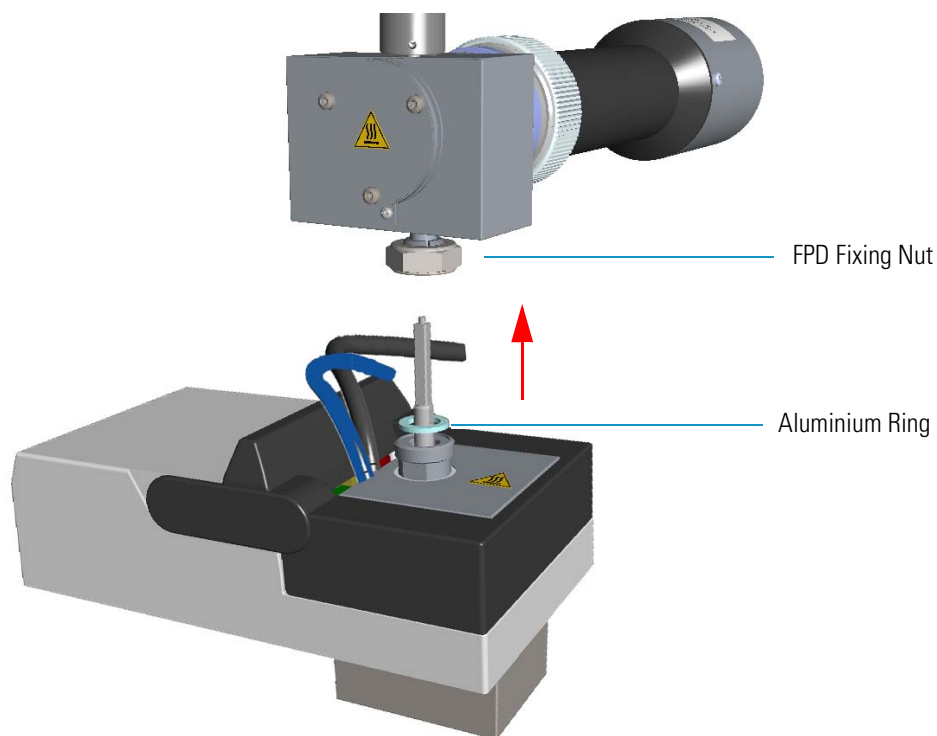

**Note** Do not lose the aluminium ring inserted between the detector head and the base body.

## Cleaning or Replacing the FPD Jet

### ❖ To clean or replace the FPD jet

1. Put the GC in standby condition.
2. Switch off the flame. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.

## 4 Performing Detectors Routine Maintenance

### Maintaining a Flame Photometric Detector (FPD)


6. Remove the FPD detector from the base.
  - a. Disconnect the signal, excitation voltage, and ignition/heating cables from the detector. See [Figure 162](#).

**Figure 162.** FPD Cable Disconnection

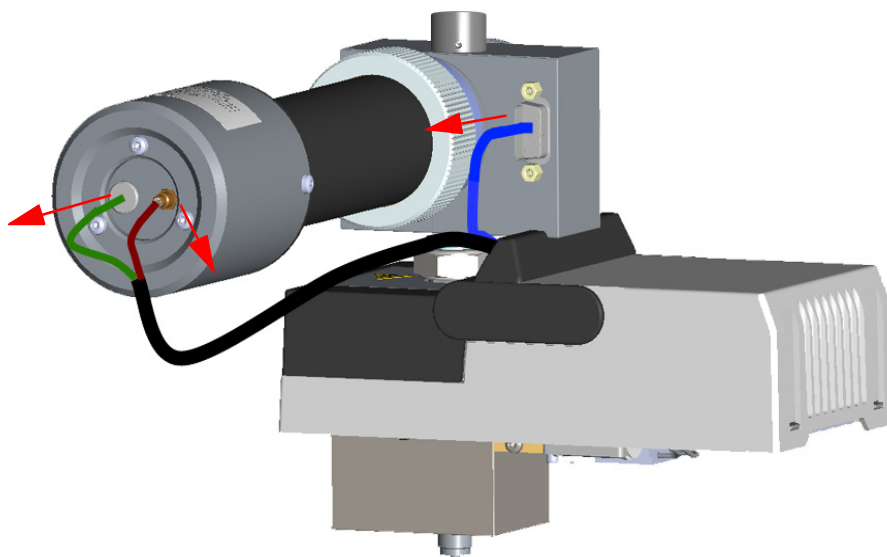

- b. Using the tool provided with the system, loosen the fixing nut on the base of the detector and remove it. See [Figure 163](#).

**Figure 163.** FPD Detector Removal

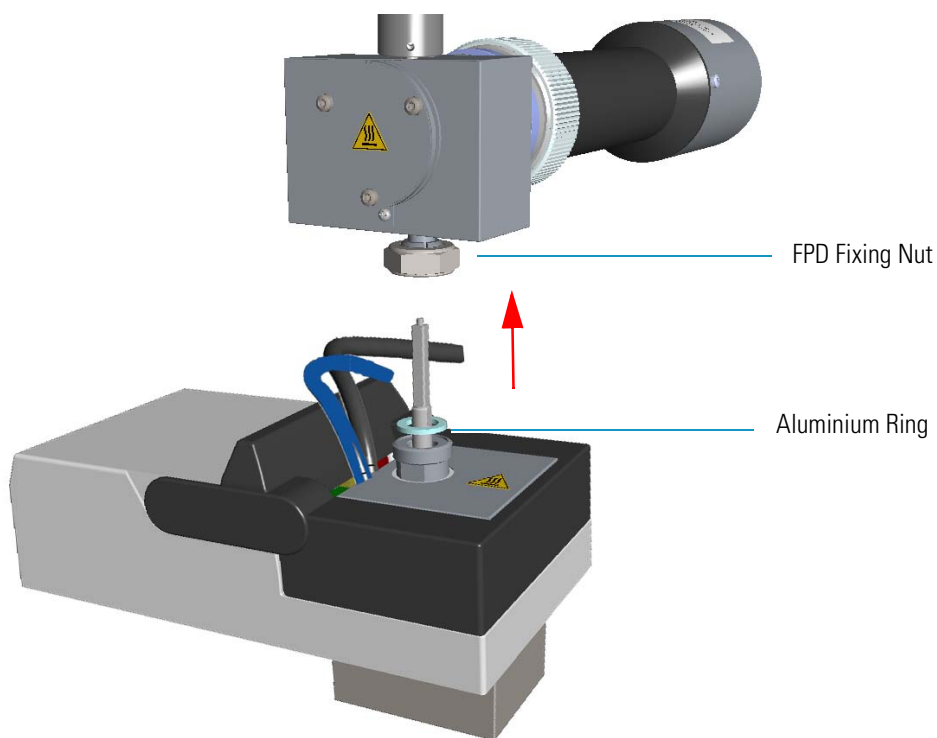


**Note** Do not lose the aluminium ring inserted between the detector head and the base body.

7. Clean the FPD Jet.


- a. Using the 5 mm wrench provided with the GC, loosen the jet and remove it from the detector base body. See [Figure 164](#).

**Figure 164.** FPD Jet Removal

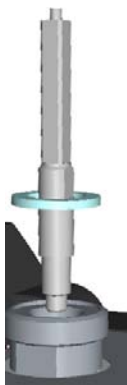

- b. Place the jet in the ultrasonic cleaner filled with liquid detergent and clean it for about five minutes.
- c. Handling the jet with forceps or tweezers, rinse the jet with distilled water, then with methanol.
- d. Place the jet on a paper towel and let it air dry. When the jet is dry, insert the jet into the detector base body and tighten it. See [Figure 165](#).

**Figure 165.** FPD Jet Reinstallation

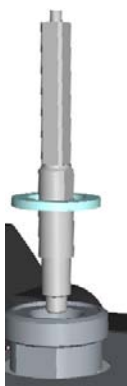

**Note** If after cleaning the physical condition of the jet does not permit its reuse, replace the jet with a new one.


8. Reinstall the FPD detector on the base.

- a. Place the detector on its base body, ensuring that the aluminium ring has been inserted in the correct position, then tighten the fixing nut. See [Figure 166](#).

## 4 Performing Detectors Routine Maintenance


### Maintaining a Flame Photometric Detector (FPD)

**Figure 166.** FPD Detector Reinstallation

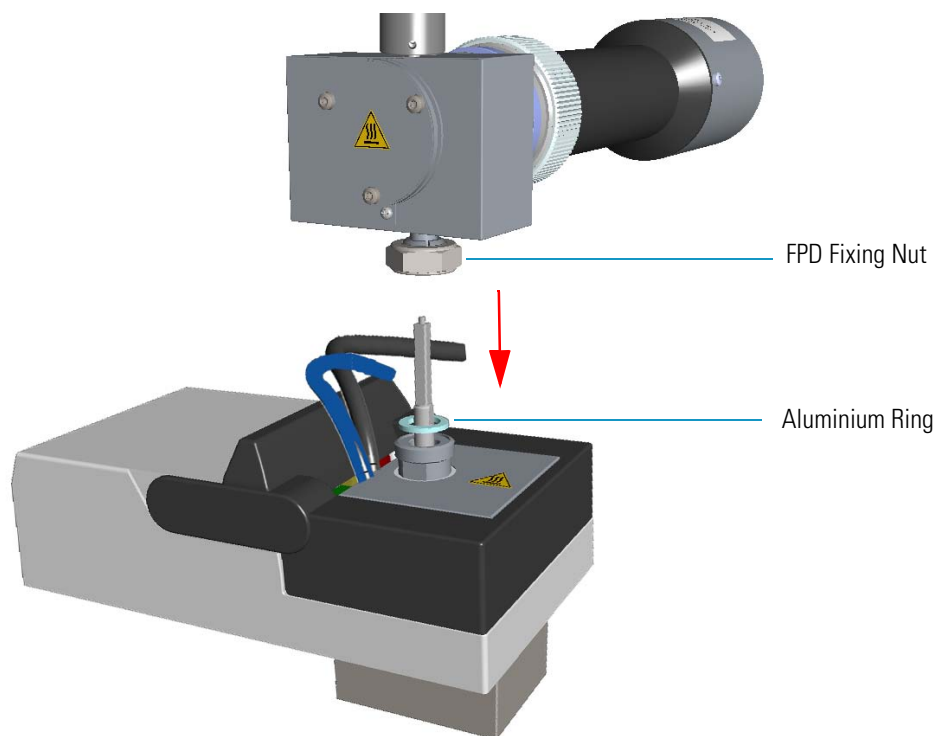

- b. Reconnect the signal, excitation voltage and ignition/heating cables to the detector. See [Figure 167](#).

**Figure 167.** FPD Cables Connection

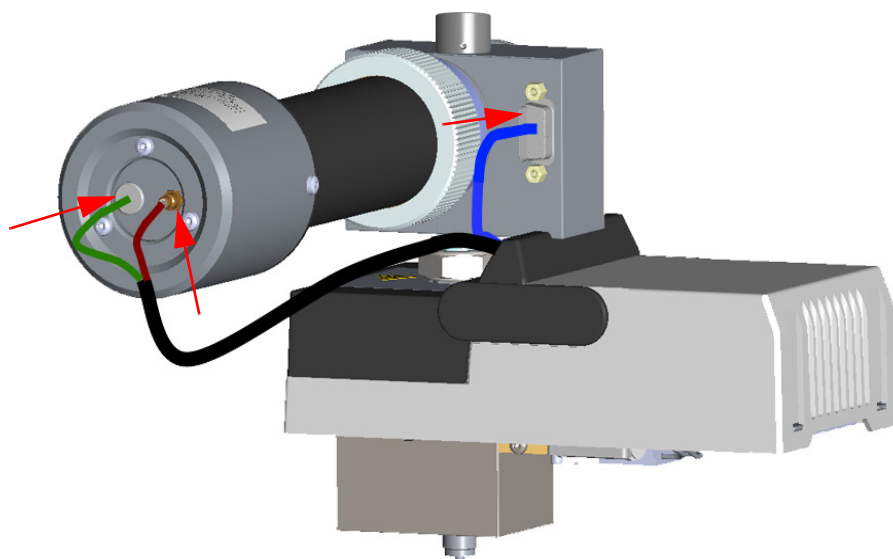

9. If other external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.
10. Power on the GC.

- a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
11. Set the normal working conditions.

## Cleaning or Replacing the FPD Interferential Filter

Before maintaining the detector, read the following precautions:



**WARNING** Carry out all the operations at low temperature to avoid burns. Therefore, before beginning maintenance, cool the detector to room temperature.



**CAUTION** When handling organic solvents you must take precautions to avoid health hazards.

### ❖ To clean or replace the FPD interferential filter


1. Put the GC in standby condition.
2. Switch off the flame. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.
6. Disconnect the signal, excitation voltage, and ignition/heating cables from the detector. See [Figure 168](#).

## 4 Performing Detectors Routine Maintenance


Maintaining a Flame Photometric Detector (FPD)

**Figure 168.** FPD Cables Disconnection

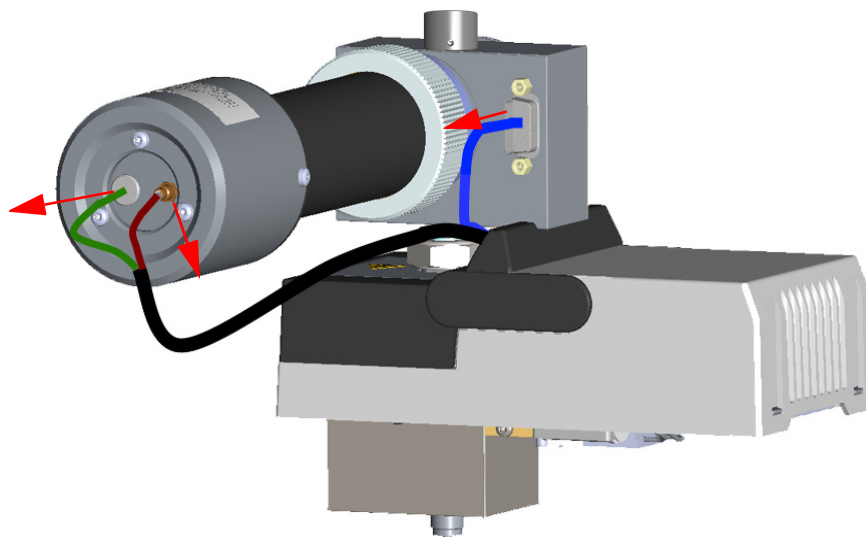

7. Remove the photomultiplier assembly and the filter. See [Figure 169](#).

**Figure 169.** Photomultiplier Assembly Removal

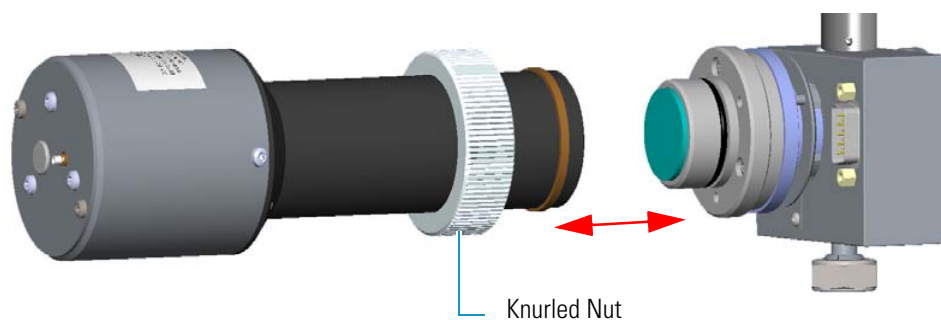


- a. Loosen the knurled nut that fixes the photomultiplier assembly and remove it from the detector body.



**CAUTION** The photomultiplier tube could damage if exposed to ambient light with the excitation voltage On. Make sure the power supply has been switched off before disconnecting the tube from the detector body.


- b. Remove the interferential filter from its housing, handling it very gently. Keep it using a clean paper towel. See [Figure 170](#).

**Figure 170.** Interferential Filter Removal

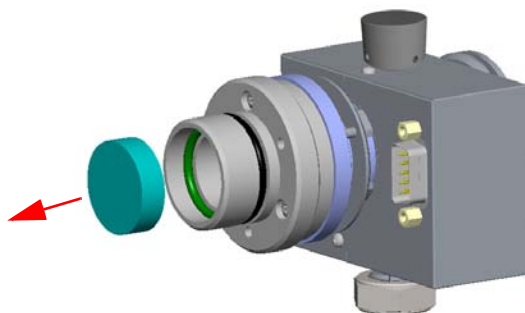

**CAUTION** Filters are fragile. Pay attention not to let the filter fall down and damage.

**CAUTION** If the filter installed is equipped with the spacer (see the kits PN 19050785 and PN 19050786), first remove the spacer, and then the interferential filter.

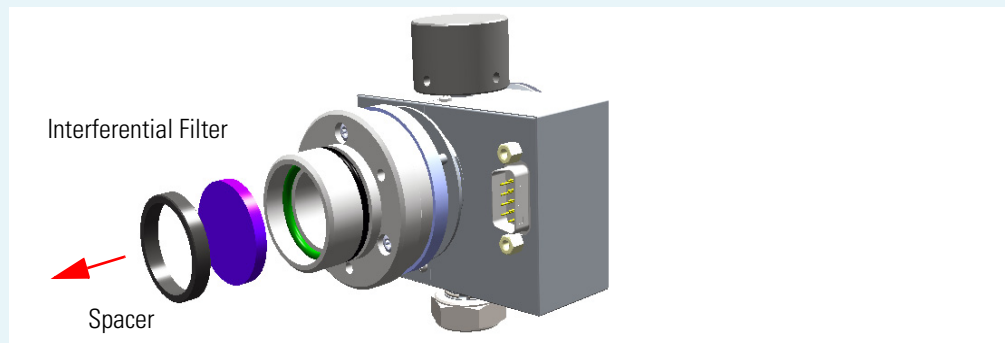

Keep them using a clean paper towel.

8. Clean the filter.
  - a. Using a clean paper towel, clean the filter on both faces.




**CAUTION** Avoid touching the filter with your fingers. If you see fingertips on the filter, clean it using GC-grade methanol and air dry before remounting.

9. Reassembly the filter and the photomultiplier assembly.
  - a. Insert the cleaned filter or a new filter into its housing. The mirror face must be oriented towards the flame. See [Figure 171](#).


## 4 Performing Detectors Routine Maintenance

### Maintaining a Flame Photometric Detector (FPD)

**Figure 171.** Filter Reinstallation

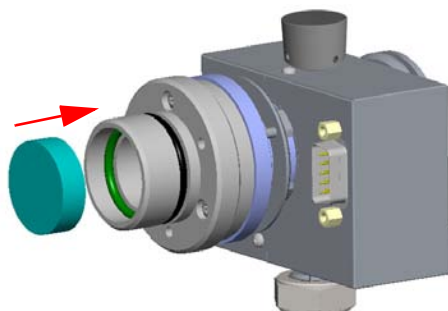


**CAUTION** If the filter installed is equipped with the spacer (see the kits PN 19050785 and PN 19050786), first insert the filter, and then the spacer.

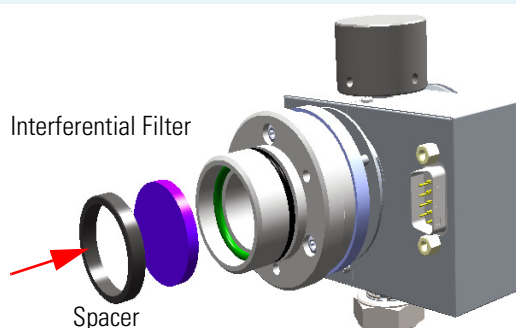

- b. Reassembly the photomultiplier assembly and the detector body, then fix them together tightening the knurled nut. See [Figure 172](#).


**Figure 172.** Photomultiplier Assembly Remounting

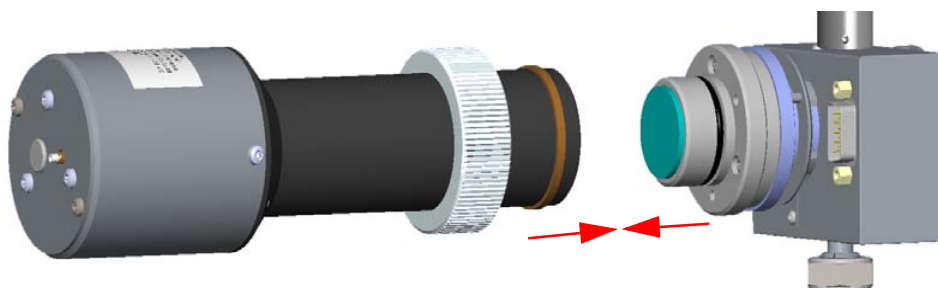

10. Reconnect the signal, excitation voltage and ignition/heating cables to the detector. See [Figure 173](#).



**Figure 173.** FPD Cables Connection

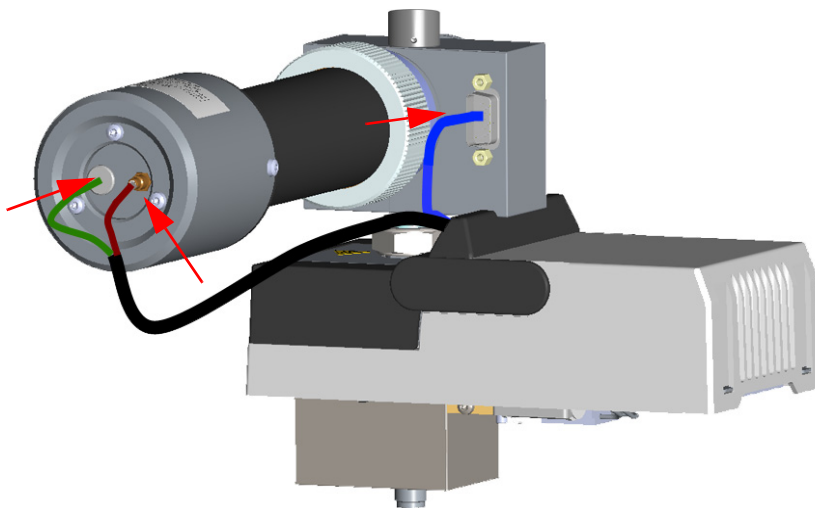

11. If other external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.
12. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
13. Set the normal working conditions.

## Replacing the FPD Ignition Glow-plug

### ❖ To replace the FPD ignition glow-plug

1. Put the GC in standby condition.
2. Switch off the flame. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.

## 4 Performing Detectors Routine Maintenance

### Maintaining a Flame Photometric Detector (FPD)


5. If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.
6. Remove the FPD detector from the base.
  - a. Disconnect the signal, excitation voltage, and ignition/heating cables from the detector. See [Figure 174](#).

**Figure 174.** FPD Cables Disconnection

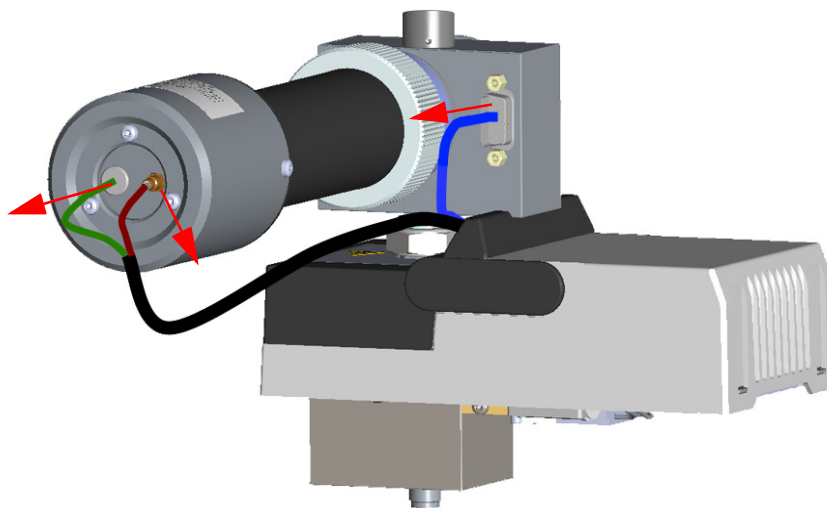

- b. Using the tool provided with the system, loosen the fixing nut on the base of the detector and remove it. See [Figure 175](#).

**Figure 175.** FPD Detector Removal

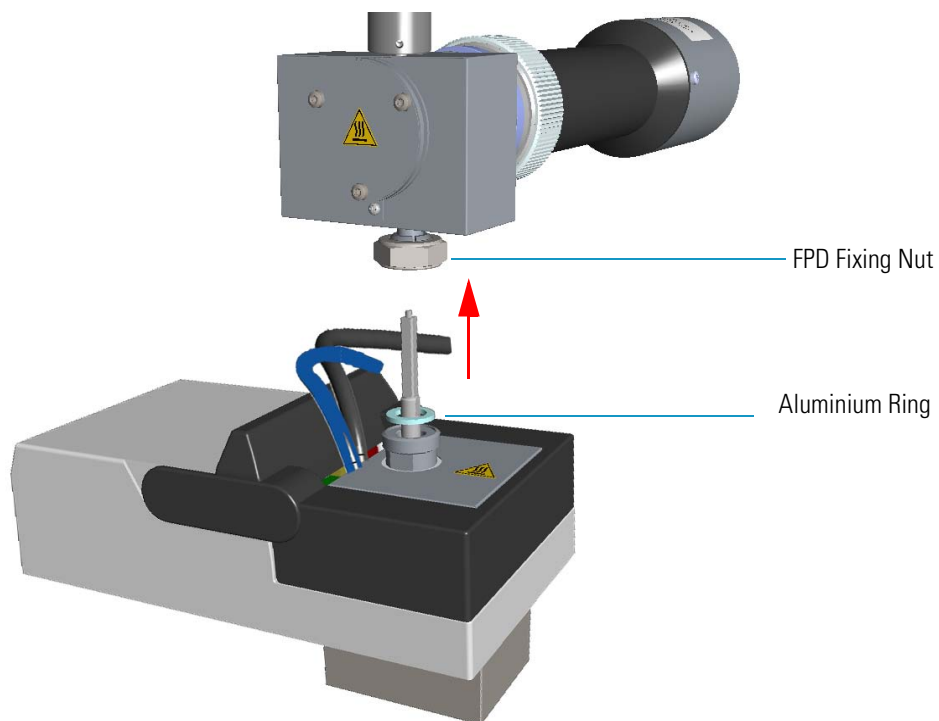

**Note** Do not lose the aluminium ring inserted between the detector head and the base body.


7. Remove the covers of the FPD detector body.



**CAUTION** Removing the covers, pay attention to the heater element and the temperature probe.


- a. Using a 5/16-in tube wrench, remove the two screws on the top of the detector body, and the front and back screws on the lower part of the detector body. See [Figure 176](#).

**Figure 176.** Detector Body Dismounting (1)

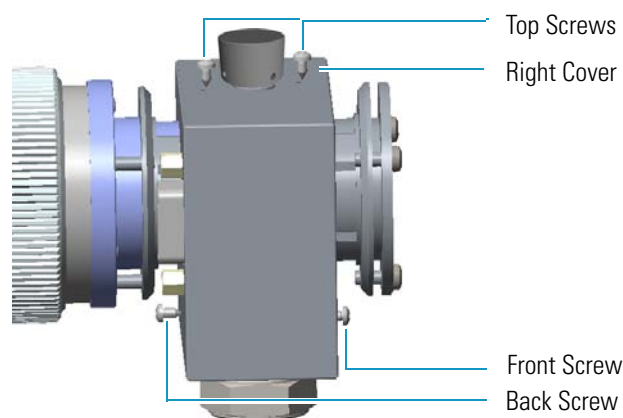

- b. Remove the left and right covers of the FPD detector body and the insulating material. The glow-plug is now accessible. See [Figure 177](#).

**Figure 177.** Detector Body Dismounting (2)

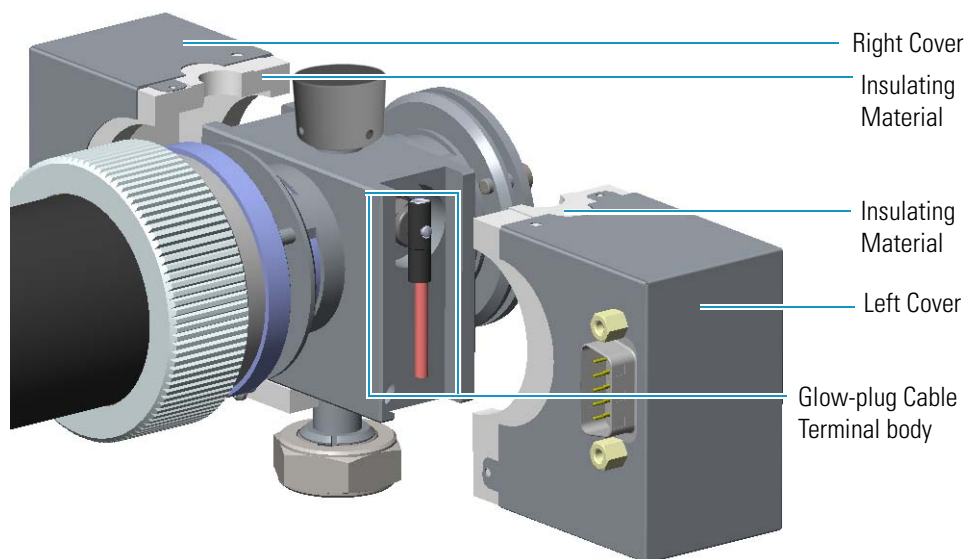


8. Replace the defective glow-plug. See [Figure 178](#).

## 4 Performing Detectors Routine Maintenance

### Maintaining a Flame Photometric Detector (FPD)

**Figure 178.** Glow-plug Replacing

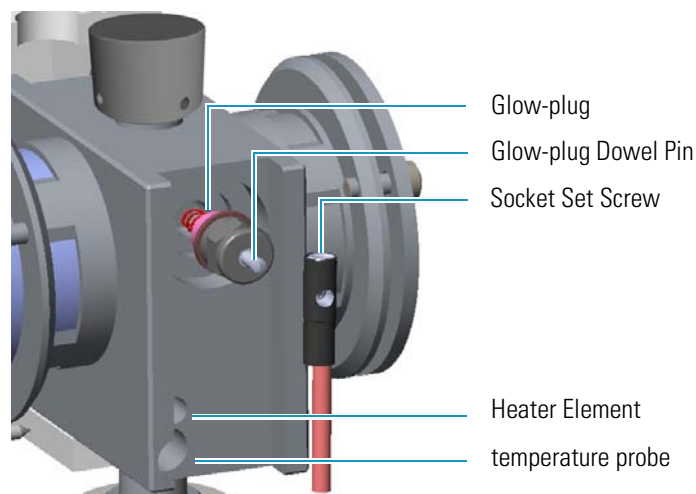


- a. Loosen the glow-plug cable socket set screw using a T6 Torxhead screwdriver. Carefully pull out the terminal body of the glow-plug cable from the dowel pin on the glow-plug element.
  - b. Using an 8-mm wrench, unscrew and remove the defective glow-plug with its washer.
  - c. Replace the glow-plug with a new one interposing the new washer.
9. Reinstall the left and right covers of the FPD detector body with the insulating material in the original position, and fix them with the four fixing screws. See [Figure 179](#).




**CAUTION** Reinstalling the covers, pay attention to the heater element and the temperature probe.

**Figure 179.** Detector Body Remounting

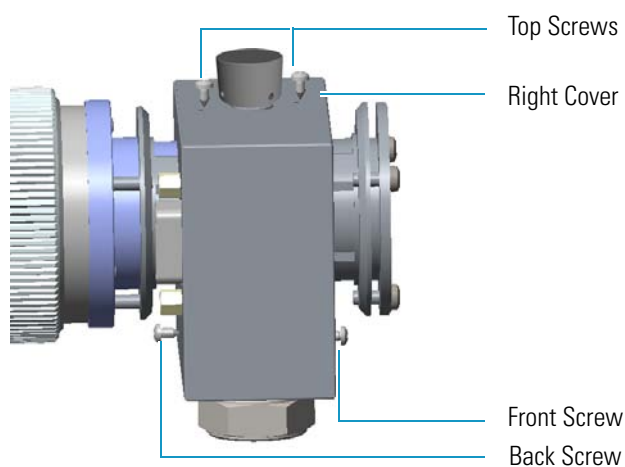


10. Reinstall the FPD detector on the base.
- a. Place the detector on its base body, ensuring that the aluminium ring has been inserted in the correct position, then tighten the fixing nut. See [Figure 180](#).

**Figure 180.** FPD Detector Reinstallation

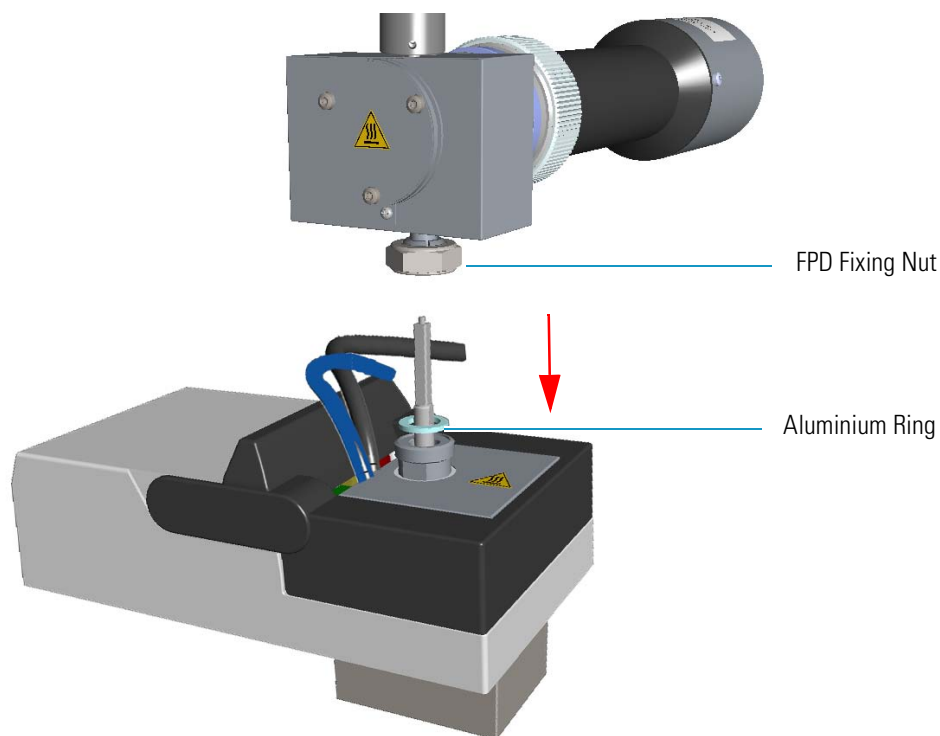

- b. Reconnect the signal, excitation voltage and ignition/heating cables to the detector.  
See [Figure 181](#).

**Figure 181.** FPD Cables Connection

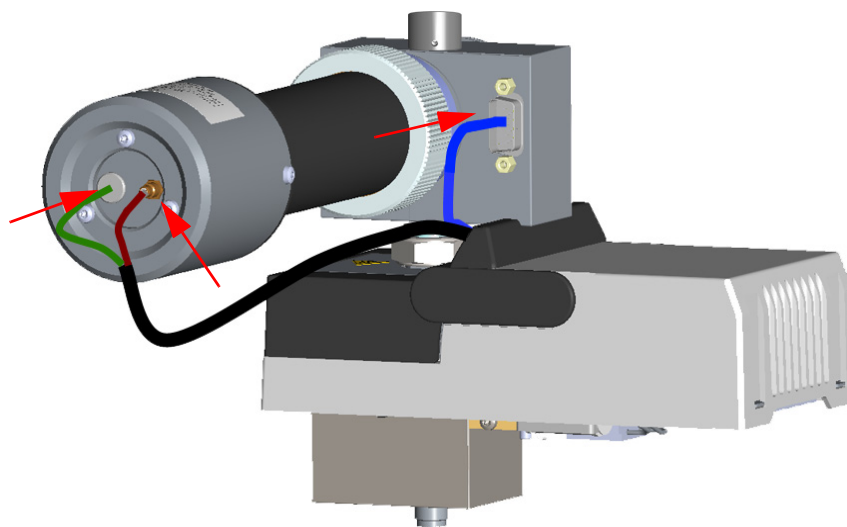

11. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
12. Power on the GC.

#### **4 Performing Detectors Routine Maintenance**

Maintaining a Flame Photometric Detector (FPD)

- a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
13. Set the normal detector working conditions.
14. Ignite the flame.

## Maintaining a Pulsed Discharge Detector (PDD)

The Pulsed Discharge Detector (PDD) does not require maintenance.



**WARNING** Under no circumstances should the detector be disassembled. The components of the detector are assembled with special tooling and held under considerable force. Disassembling of the detector may present a safety hazard and will result in its destruction.





## GC Main Frame Advanced Maintenance

This chapter describes TRACE 1300/TRACE 1310 GC components that do not require routine maintenance, but they need to be removed or replaced.

### Contents

- [Removing/Replacing the GC Top Cover](#)
- [Removing/Replacing the GC Left Side Panel](#)
- [Removing/Replacing the GC Right Side Panel](#)
- [Removing/Replacing the GC Back Cover](#)
- [Removing/Replacing the GC Front Door Cover](#)
- [Removing/Replacing the Electronic Module](#)
- [Replacing the Oven Heater Baffle](#)
- [Replacing the Oven Heater Temperature Sensor](#)
- [Replacing the Oven Motor](#)
- [Replacing the Flap Motor](#)

## Removing/Replacing the GC Top Cover

Perform this operation for accessing the internal compartments of the GC.



**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.


❖ **To remove/replace the top cover accessing the top parts of the GC**

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Power off the GC.
  - a. Close the gas supplies.
  - b. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - c. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
5. Remove the autosampler if present.
6. Remove the cover.
  - a. Using a T20 Torxhead screwdriver, unscrew the four screws that held it in place. See [Figure 182](#) and [Figure 183](#).

**Figure 182.** Top Cover Removal (1)

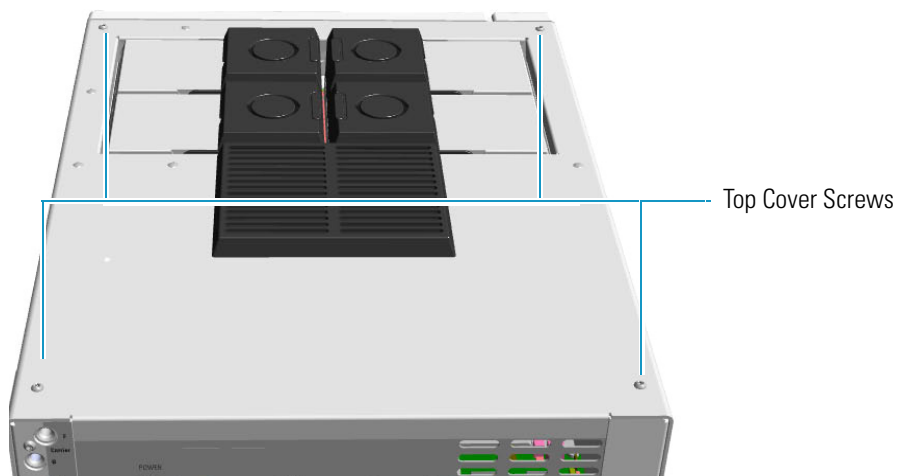

- b. Lift the top panel up and off the GC.

**Figure 183.** Top Cover Removal (2)



7. Replace the top cover proceeding in the reverse order in which it was removed.

## Removing/Replacing the GC Left Side Panel

❖ **To remove the left side panel with the side panel molded of the GC**

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

## 5 GC Main Frame Advanced Maintenance

### Removing/Replacing the GC Left Side Panel


3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
5. Remove the panel.
  - a. Open the front door of the GC. Using a T20 Torxhead screwdriver, unscrew the left side panel screw from the interior front panel. See [Figure 184](#).

**Figure 184.** Left Side Panel Fixing Screw

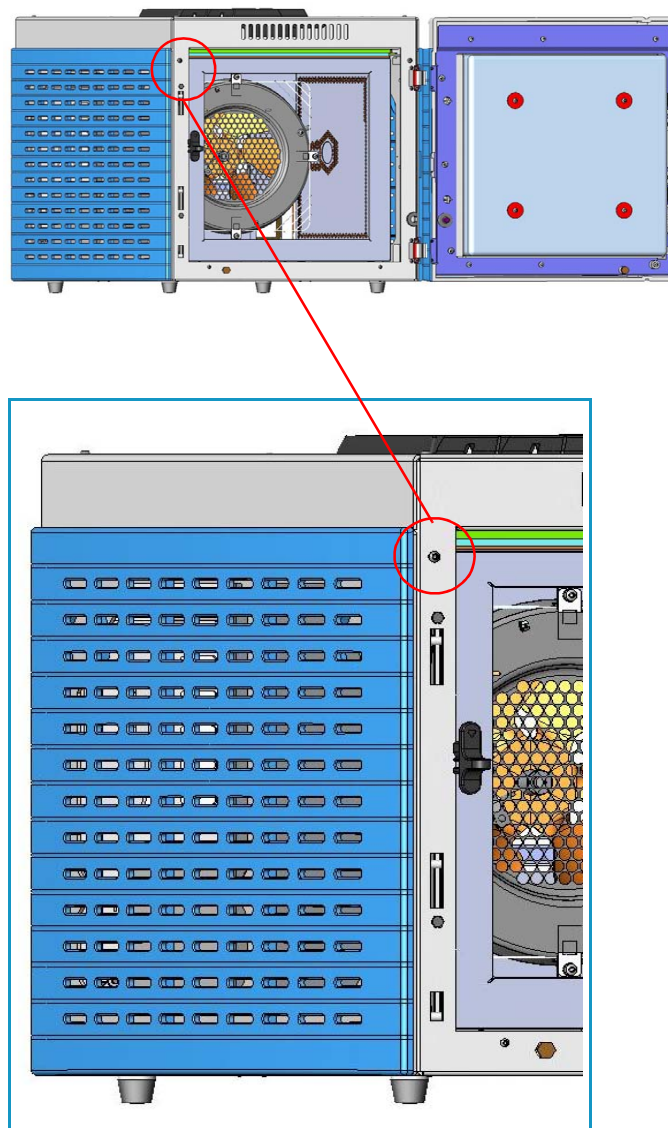

- b. Slide the panel towards the back of the instrument up to the stop.
- c. Remove the panel pulling it outwards being aware that the ground wire is attached to the panel. See [Figure 185](#).

**Figure 185.** Left Panel Removal

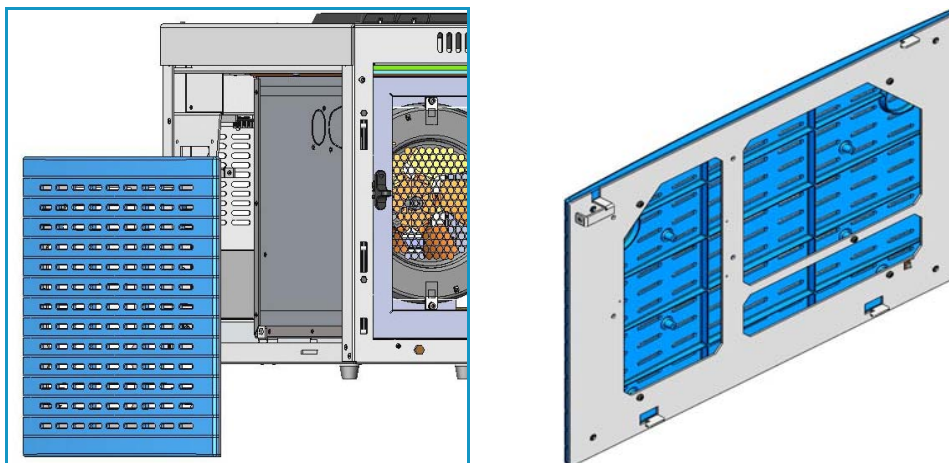

- d. Unplug the ground wire from the panel.

**Note** Pay attention to the positioning of the ground wire plug, so it can be reconnected in the same way it was removed.

6. Replace the panel.
  - a. Plug the ground wire to the panel.
  - b. Reinstall the panel proceeding in the reverse order in which it was removed.

## Removing/Replacing the GC Right Side Panel

### ❖ To remove the right side panel with the side panel molded of the GC

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.

## 5 GC Main Frame Advanced Maintenance

### Removing/Replacing the GC Right Side Panel


4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
5. Remove the AI/AS autosampler, if present.
6. Remove the panel.
  - a. Open the front door of the GC. Using a T20 Torxhead screwdriver, unscrew the right side panel screw from the interior front panel. See [Figure 186](#).

**Figure 186.** Right Panel Fixing Screw

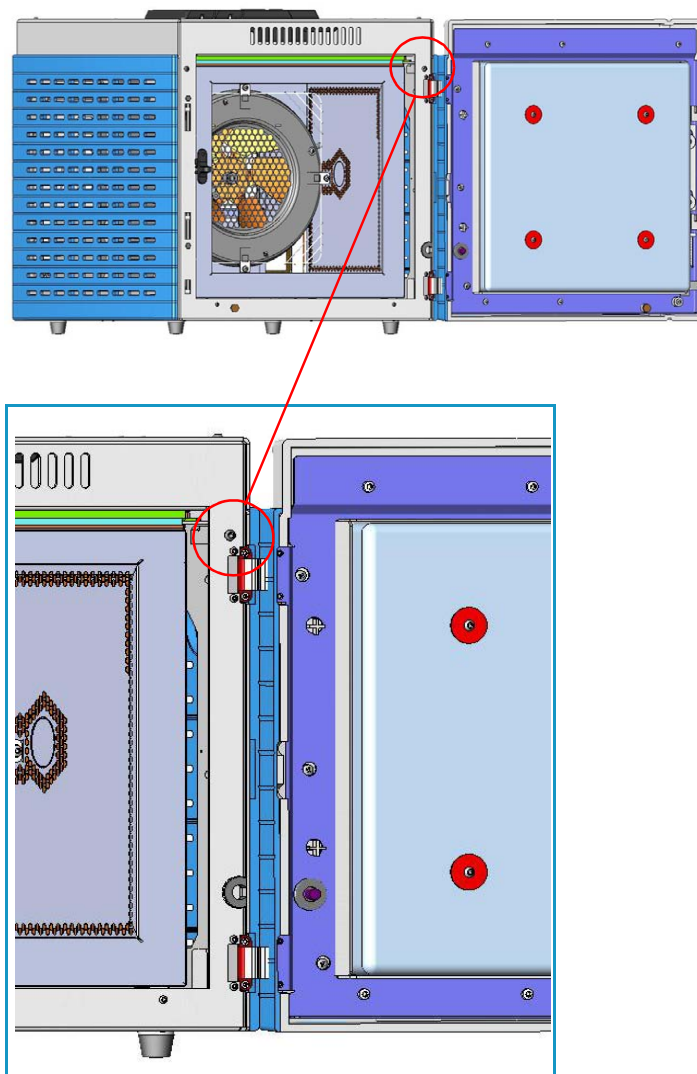

- b. Slide the panel towards the back of the instrument up to the stop.
  - c. Remove the panel pulling it outwards being aware that the ground wire is attached to the panel. See [Figure 187](#).

**Figure 187.** Right Panel Removal

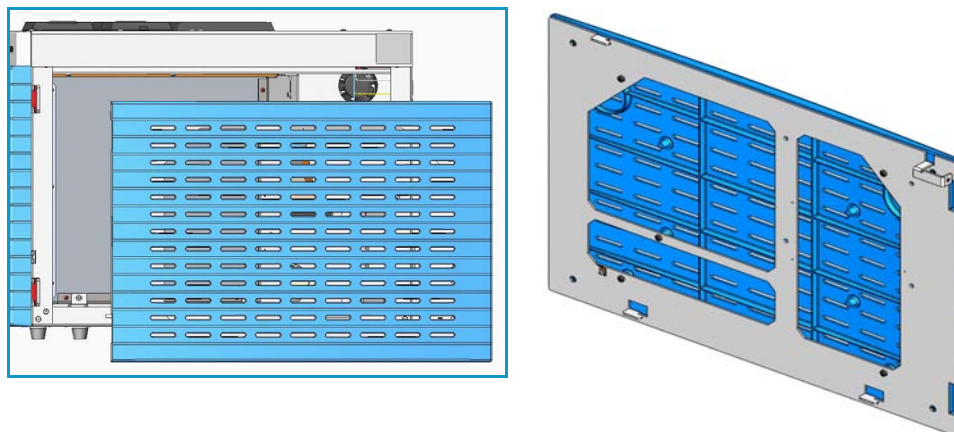

- d. Unplug the ground wire from the panel.

**Note** Pay attention to the positioning of the ground wire plug, so it can be reconnected in the same way it was removed.

7. Replace the panel.
  - a. Plug the ground wire to the panel.
  - b. Reinstall the panel proceeding in the reverse order in which it was removed.

## Removing/Replacing the GC Back Cover

Remove the back cover for accessing the compartment on the back of the GC.



**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.

### ❖ To remove the back cover

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.




## 5 GC Main Frame Advanced Maintenance

### Removing/Replacing the GC Back Cover


4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
5. Remove the cover.
  - a. Using a T20 Torxhead screwdriver, unscrew the four screws that secure the back cover to the GC. See [Figure 188](#).

**Figure 188.** Back Cover Removal

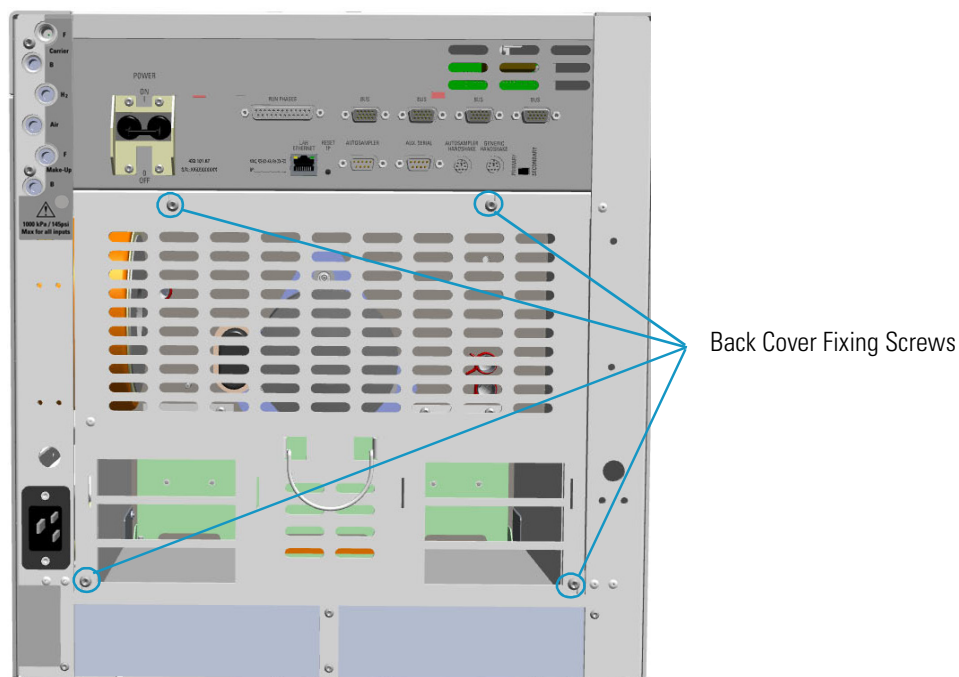

- b. Lift the cover off using the cover handle. Be aware that the ground wire is attached to the back cover. See [Figure 189](#).

**Figure 189.** Back Cover

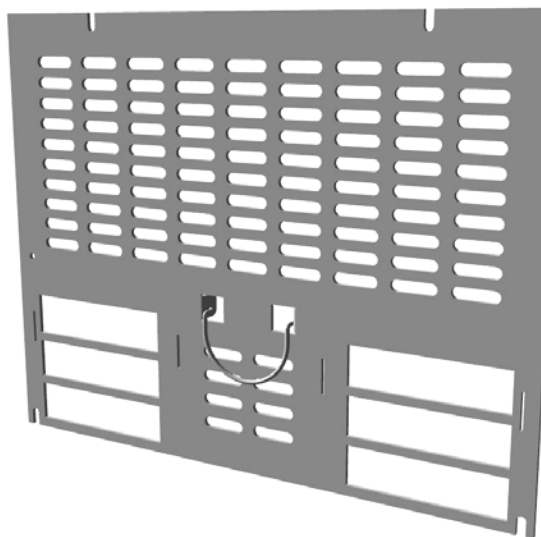



**Note** Pay attention to the positioning of the ground wire plug, so it can be reconnected in the same way it was removed.

6. Reconnect the ground wire to the back cover terminal.
7. Replace the cover proceeding in the reverse order in which it was removed.

## Removing/Replacing the GC Front Door Cover

This section provides the instruction for removing or replacing the cover of the front door.



**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.

### ❖ To remove the cover of the front door of the GC

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
5. Remove the door cover.
  - a. Open the front door and look for the fixing screws that secure the cover and the door handle to the door. See [Figure 190](#).

## 5 GC Main Frame Advanced Maintenance

Removing/Replacing the GC Front Door Cover

**Figure 190.** Front Door Back View

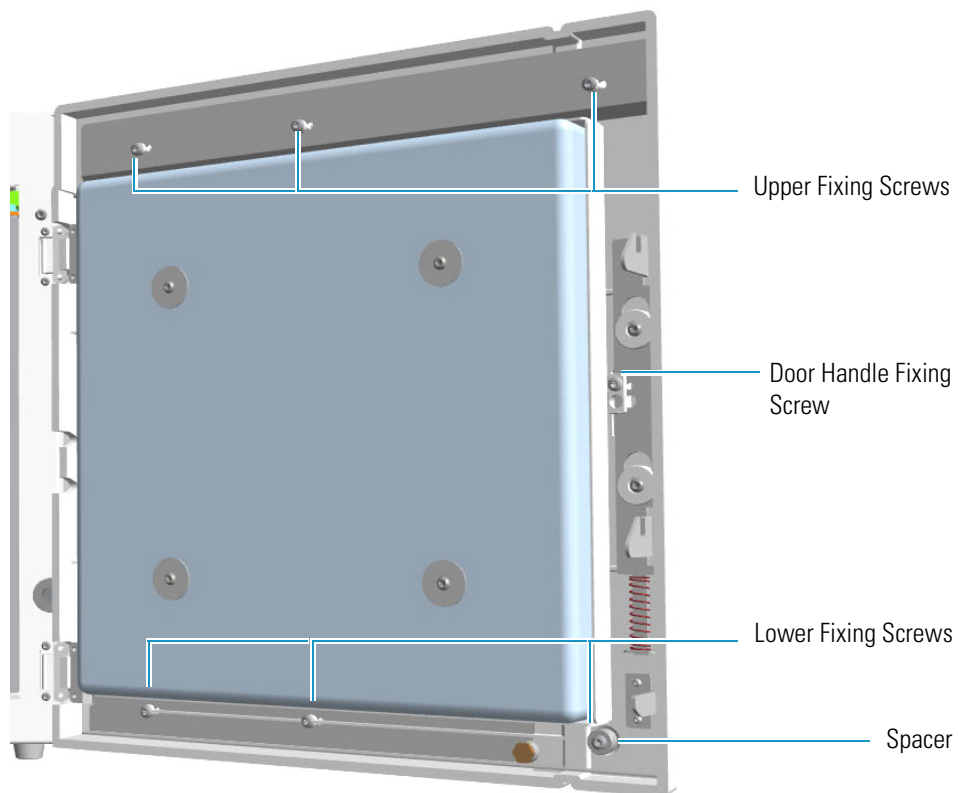


- b. Using a T20 Torxhead screwdriver remove the screw that secure the door handle to the front door.
- c. Pull the door handle out from the front door. See [Figure 191](#).



**IMPORTANT** Save the door handle because it must be reused.


**Figure 191.** Door Handle Removal

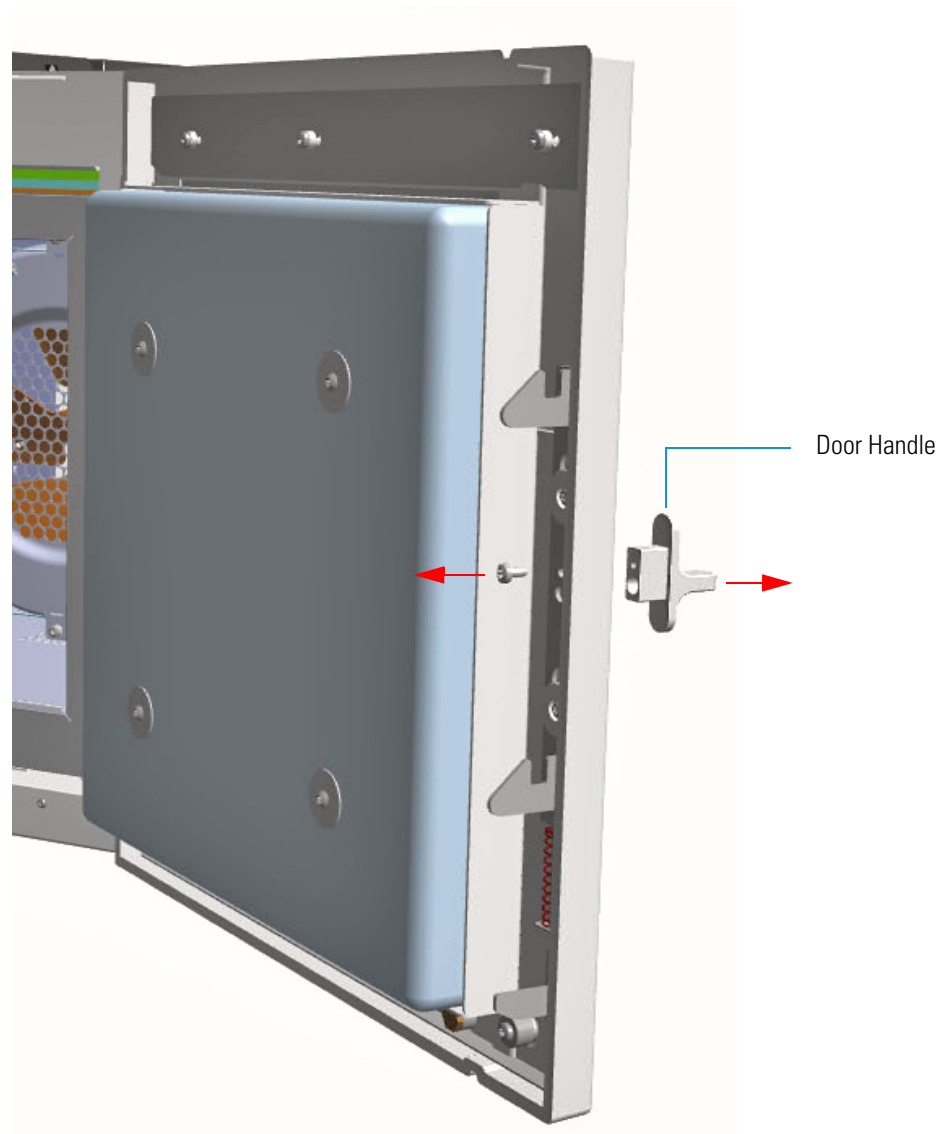

- d. Using a T20 Torxhead screwdriver remove the three upper and the three lower screws that secure the cover to the front door. See [Figure 192](#).


**Note** The lower screw on the lower right corner is screwed into a spacer.

**Figure 192.** Door Cover Fixing Screws Removal

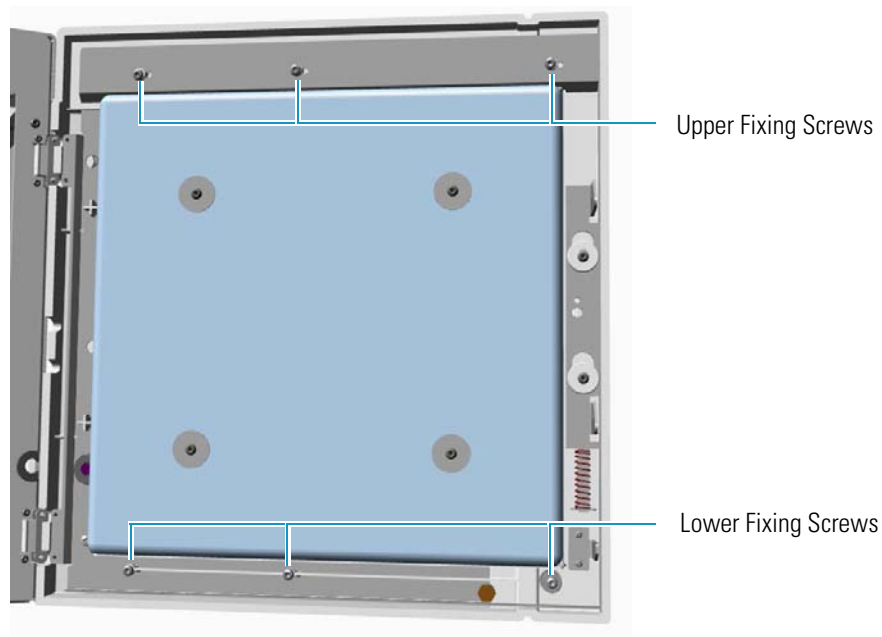

- e. Carefully pull the door cover off. Pay attention to the cables that connect the status panel to the internal section of the door. See [Figure 193](#).


**Figure 193.** Front Door Removal

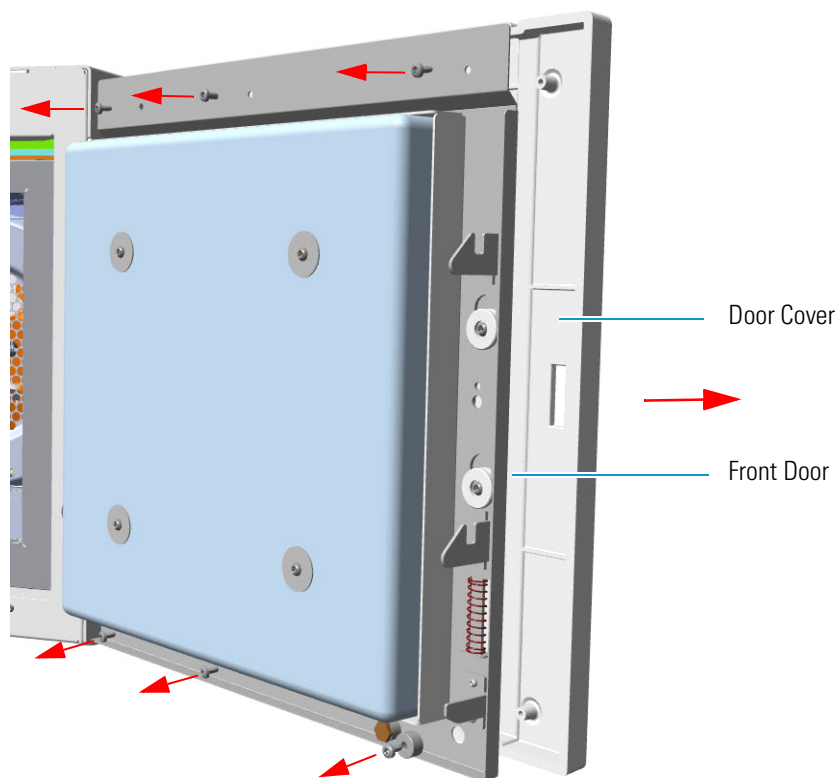


- f. Disconnect the cables from the connector located on the front of the door.

- See [Figure 194](#) in case of TRACE 1300.
- See [Figure 195](#) in case of TRACE 1310.

**Figure 194.** TRACE 1300 Front Door Cables Removal

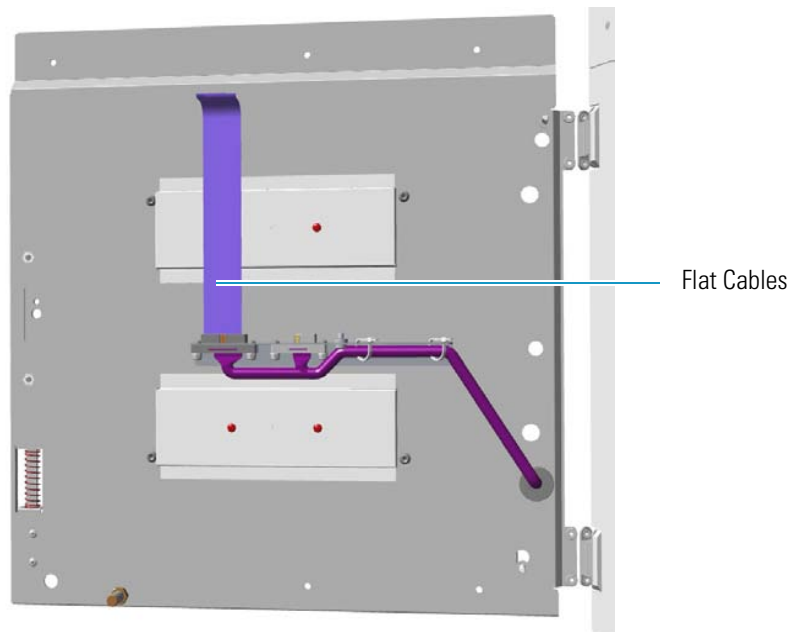

**Figure 195.** TRACE 1310 Front Door Cables Removal

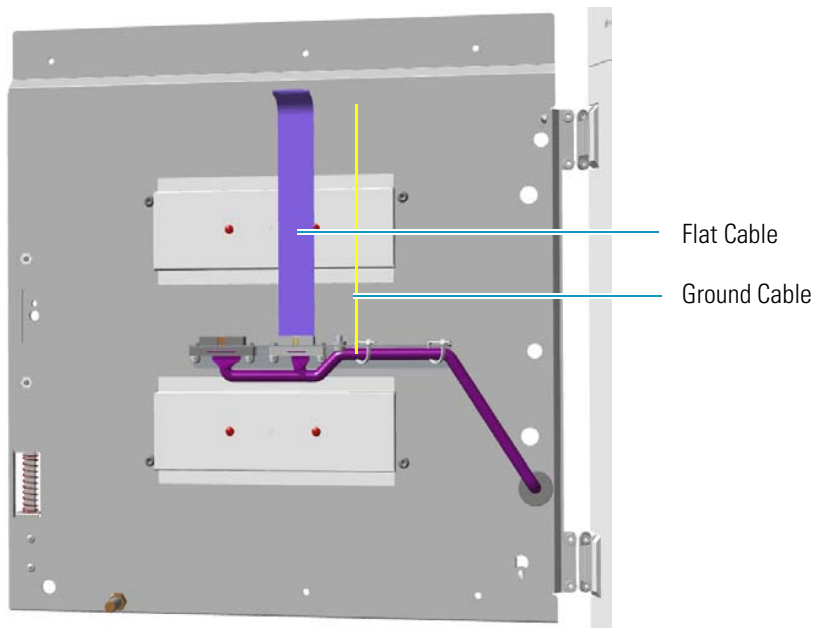

6. Replace the cover proceeding in the reverse order in which it was removed.
7. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.

8. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
9. Set the normal injector, detector, and GC working conditions.

## Removing/Replacing the Electronic Module

The Electronic Module contains power and electronic circuits for the control of the instrument. A proper module is installed according to 120 Vac or 230 Vac power supply.




**WARNING** This operation must be carried out only by authorized and trained Thermo Fisher Scientific technical personnel.

### ❖ To remove/replace the electronic module

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.
3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
5. Unplug the cables connected to the external interface of the Electronic Module.
6. Remove the autosampler if present.
7. Remove the top cover.
  - a. Use a T20 Torxhead screwdriver to remove the screws at the sides of the top cover. See [Figure 196](#).

**Figure 196.** Top Cover Removal (1)

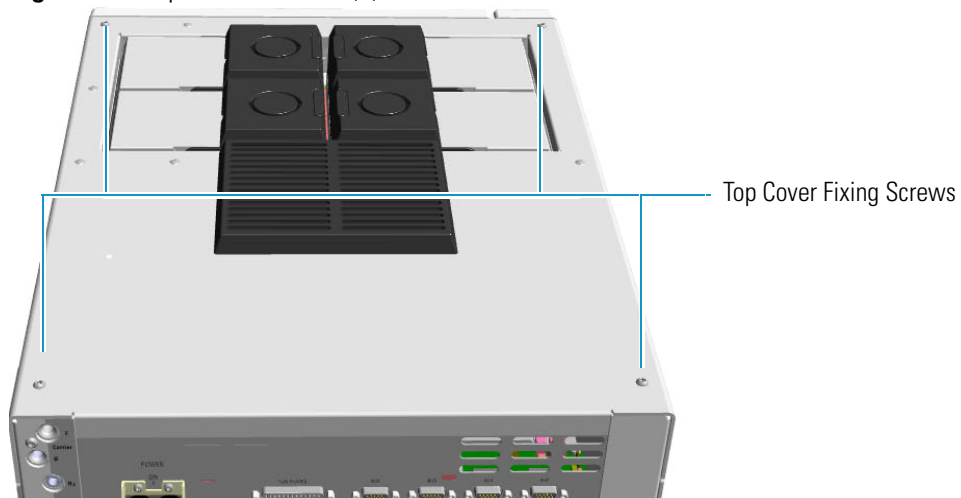

- b. Slide the top panel toward the back of the instrument and lift it off. See [Figure 197](#).


**Figure 197.** Top Cover Removal (2)

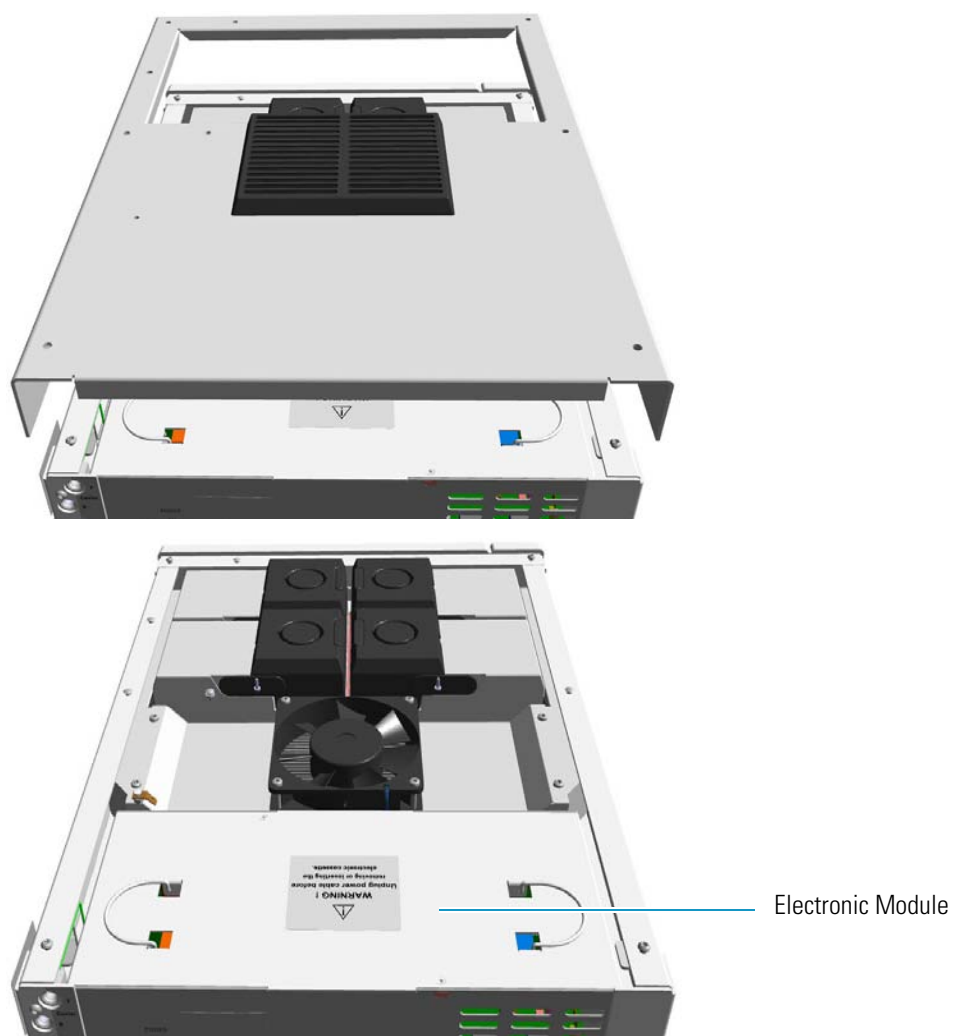

8. Remove the Electronic Module.

- a. Using the two handles on the top cover, pull out module from its housing. See [Figure 198](#).


**Figure 198.** Electronic Module Extraction

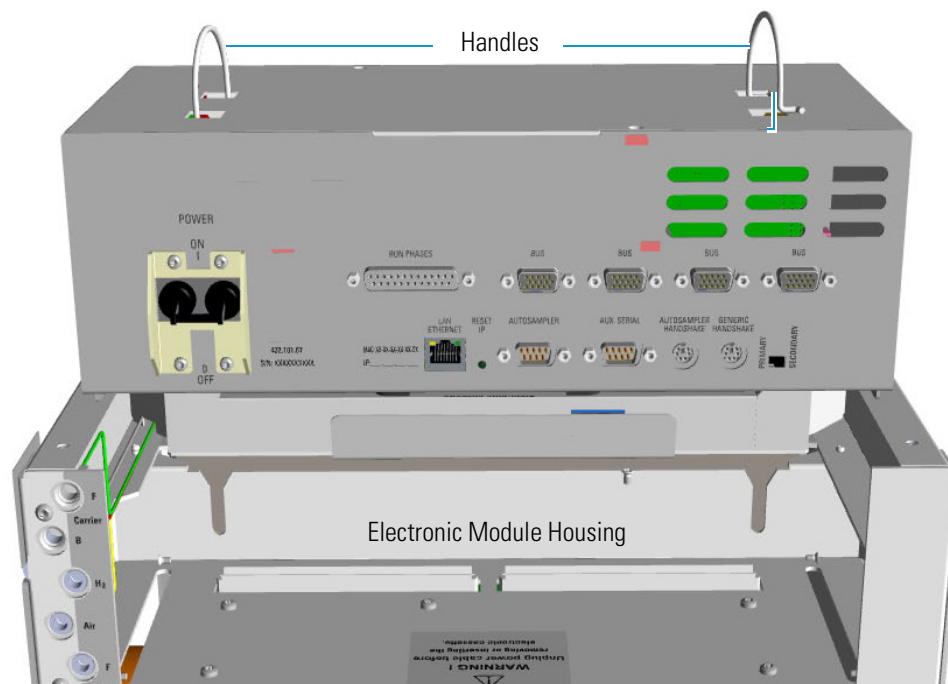

9. Replace the Electronic Module.

- a. Guide the electronic module into its housing.
- b. Push down the module. Be sure to plug the terminal contacts of the mother board into the two slots of the backplane board. See [Figure 199](#) and [Figure 200](#).

**Figure 199.** Electronic Module Replacement (1)

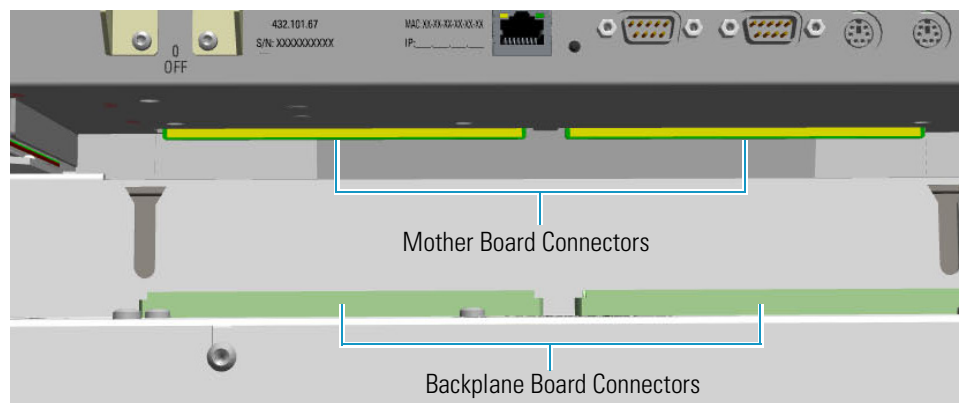






❖ **To replace the oven heater baffle**

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
5. Remove the back cover.
  - a. Using a T20 Torxhead screwdriver, remove the four screws that secure the back cover to the GC. See [Figure 201](#).

**Figure 201.** Back Cover Removal

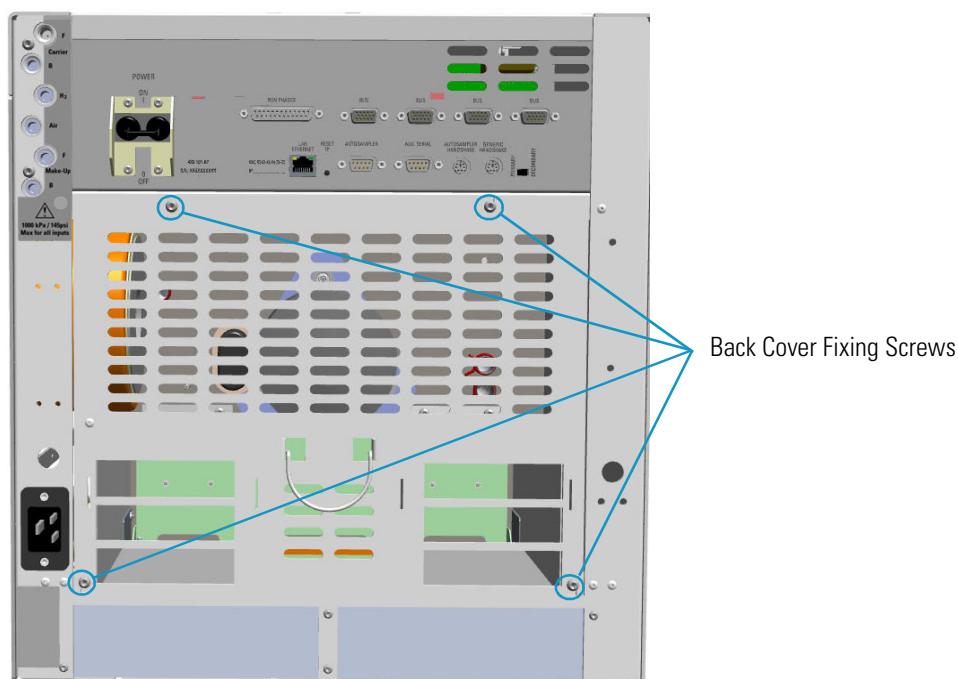


- b. Lift the cover off using the cover handle. Be aware that the ground wire is attached to the back panel.

**Note** Pay attention to the positioning of the ground wire plug, so it can be reconnected in the same way it was removed.

6. Remove the analytical column.
  - a. Open the front door of the GC.
  - b. Loosen the retaining nut from the injector and detector fitting on the upper interior wall of the GC oven.
  - c. Remove the analytical column with its nut and ferrule from the bottom of the injector and the detector.
7. Remove the oven heater baffle.

The oven heater baffle is attached to the oven wall by four adjustable tabs. These tabs are attached to the oven heater baffle with four Torx screws. See [Figure 202](#).

**Figure 202.** Oven Heater Baffle

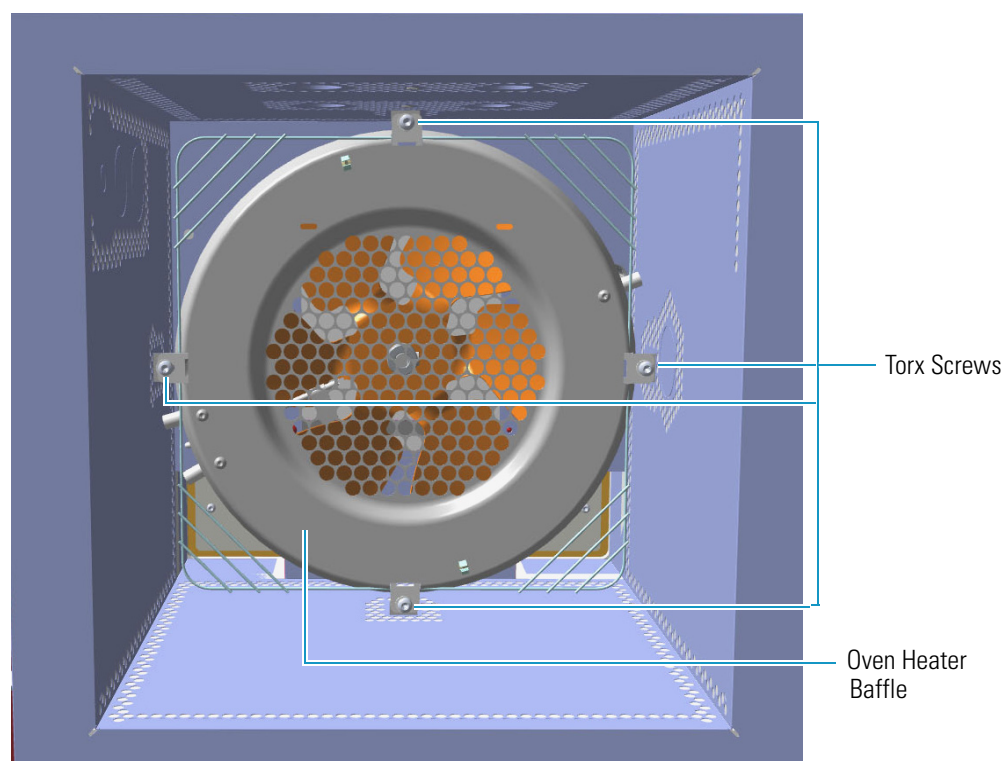


The heater and temperature sensor wires are connected to the backplane board passing through the three pipes provided on the oven back wall.

- The temperature sensor wire is plugged into the connector marked **J4** (Oven PT100)
- The heater wires are plugged into the connector **J23**, **J24**, and **J25** (Oven Heater).

The pipes are plugged by insulation material. See [Figure 203](#).

**Figure 203.** Pipes

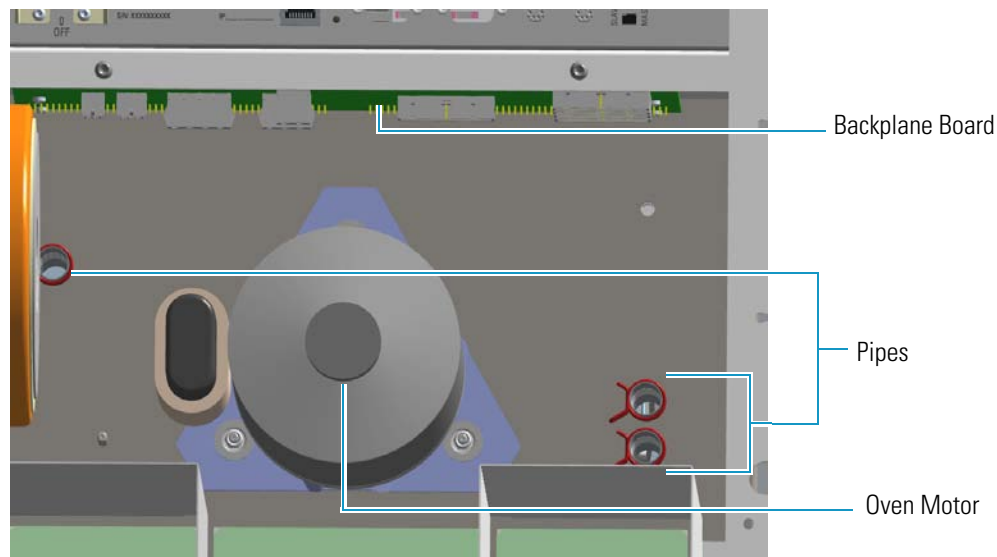

- a. Unplug the heater and temperature sensor wires from the relevant connector on the backplane board.
  - b. Loosen and remove the four Torx screws that secure the tabs to the oven heater baffle. See [Figure 202](#).
8. Pull out the oven heater baffle from the oven, paying attention to guide the heater and temperature sensor wires through the three pipes on the back oven wall.



**IMPORTANT** During this operation, parts of the insulation material drops down into the oven. This material is reused later when a new Oven Heater Baffle is installed.

9. Replace the oven heater baffle in the reverse order in which it was removed.
  - a. Guide the heater and temperature sensor wires through the pipes provided on the oven back wall.
  - b. Reconnect the temperature sensor and heater wires to the relevant connector **J4**, **J23**, **J24**, and **J25** on the backplane board.
  - c. Plug the pipes with the insulation material.
  - d. Secure the oven heater baffle to the oven wall with the four Torx screws in the tabs.

The tab adjustment screws should be loose. Carefully center the oven heater baffle inside of the oven. Be sure that the center hole in the oven heater baffle screen is centered over the blower motor shaft. The tolerance is very small.

The oven heater baffle will not move left or right, nor up or down more than a couple of millimeters. However, be sure that it is centered in the oven over the blower wheel.

This will ensure that whenever the GC is turned on, the blower wheel will not be touching the heater element mounted to the oven heater baffle.

10. Reinstall the back cover.
  - a. Reconnect the ground wire to the back cover terminal.
  - b. Replace the cover proceeding in the reverse order in which it was removed.
11. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
12. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
13. Close the front door of the GC. Be sure that the oven blower is running.
  - a. Verify that the heater heats, and the fan rotates without contacting the heater element.
14. Open the front door of the GC.
15. Reinstall the column.
16. Set the normal injector, detector and GC working conditions.

## Replacing the Oven Heater Temperature Sensor

Removing the oven heater temperature sensor requires removing the complete oven heater baffle.



**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.

### ❖ To replace the oven heater temperature sensor


1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.

- c. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
4. Remove the back cover.
  - a. Using a T20 Torxhead screwdriver remove the four screws that secure the back cover to the GC. See [Figure 204](#).

### Figure 204. Back Cover Removal

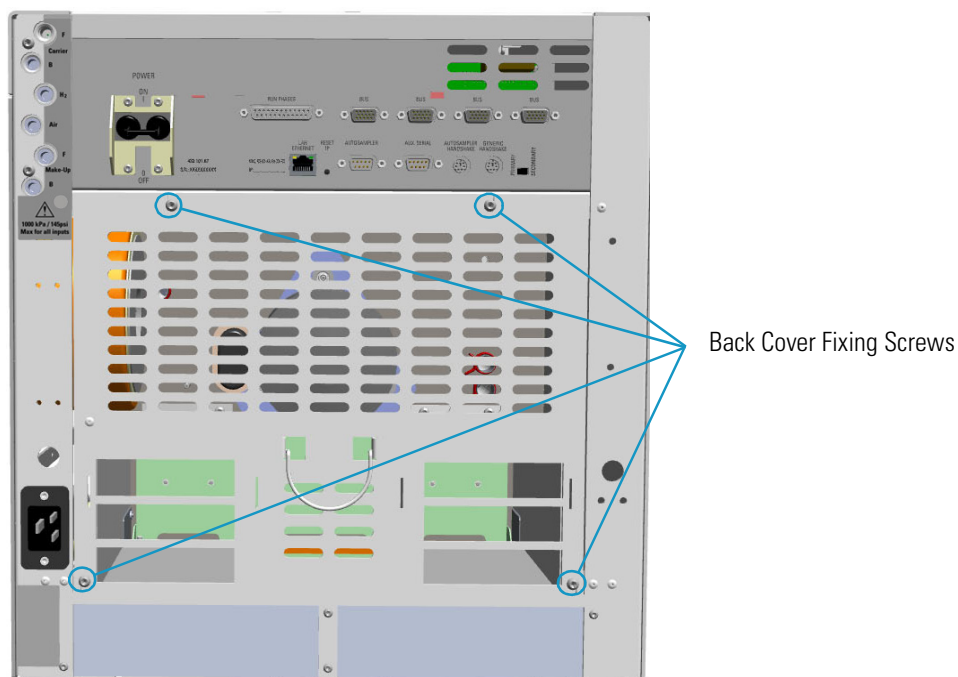


- b. Lift the cover off using the cover handle. Be aware that the ground wire is attached to the back panel.

**Note** Pay attention to the positioning of the ground wire plug, so it can be reconnected in the same way it was removed.

5. Remove the analytical column.
  - a. Open the front door of the GC.
  - b. Loosen the retaining nut from the injector and detector fitting on the upper interior wall of the GC oven.
  - c. Remove the analytical column with its nut and ferrule from the bottom of the injector and the detector.
6. Remove the oven heater baffle.

The oven heater baffle is attached to the oven wall by four adjustable tabs. These tabs are attached to the oven heater baffle with four Torx screws. See [Figure 205](#).

**Figure 205.** Oven Heater Baffle

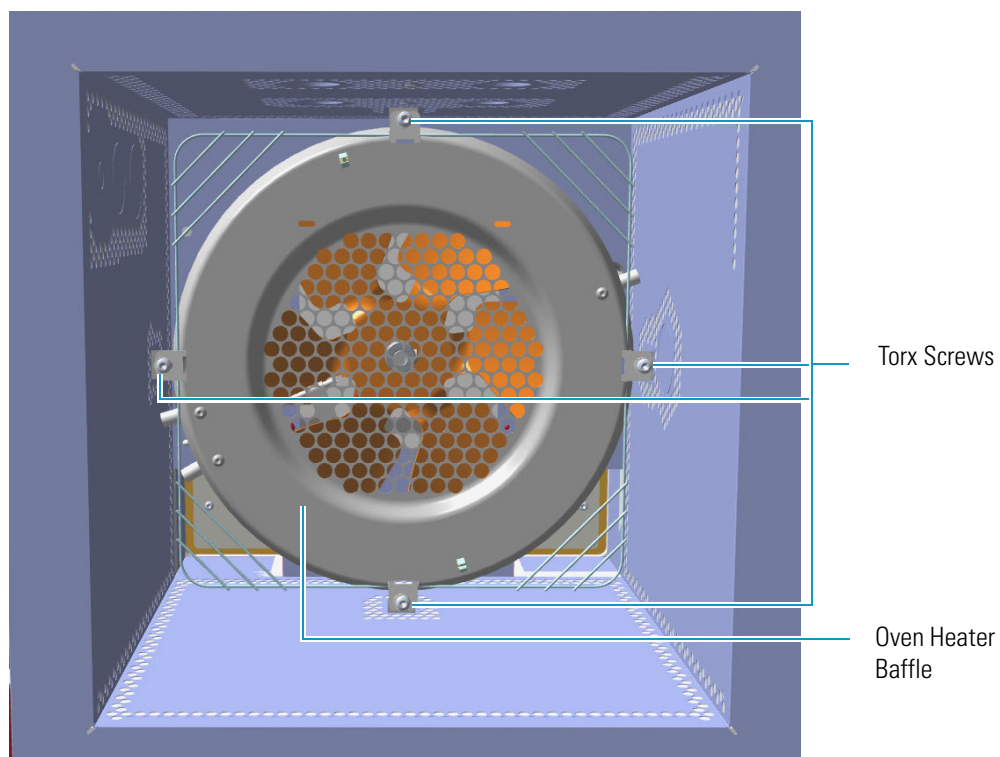


The heater and temperature sensor wires are connected to the backplane board passing through the three pipes provided on the oven back wall.

- The temperature sensor wire is plugged into the connector marked **J4** (Oven PT100)
- The heater wires are plugged into the connector **J23**, **J24**, and **J25** (Oven Heater).

The pipes are plugged by insulation material. See [Figure 206](#).

**Figure 206.** Pipes

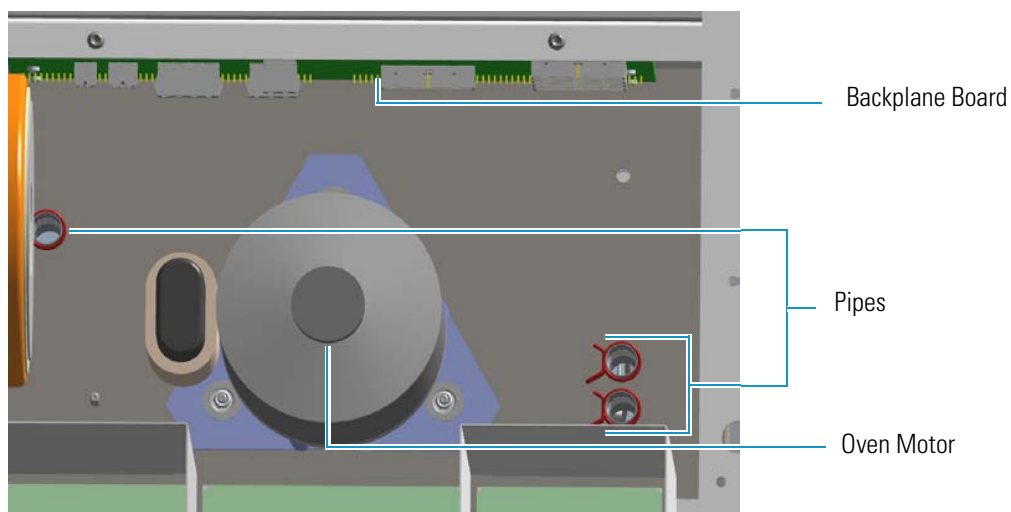




- a. Unplug the heater and temperature sensor wires from the relevant connector on the backplane board.
- b. Loosen and remove the four Torx screws that secure the tabs to the oven heater baffle. See [Figure 205](#) on [page 217](#).
- c. Place the heater baffle on the roof of the oven.



**IMPORTANT** During this operation, parts of the insulation material drops down into the oven. This material is reused later when a new Oven Heater Baffle is installed.

7. Remove the oven heater temperature sensor. See [Figure 207](#).

**Figure 207.** Oven Heater Temperature Sensor

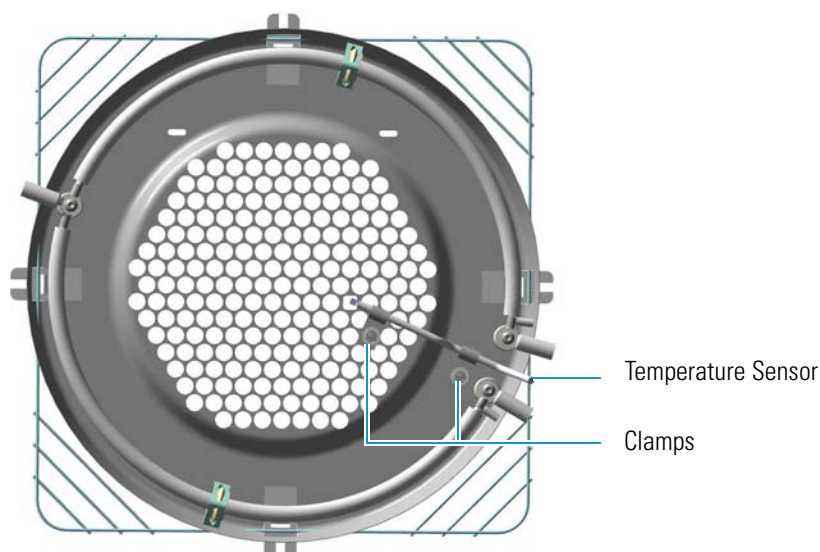

- a. Loosen the two clamps that attaches the oven heater temperature sensor to the oven heater baffle.
- b. Remove the temperature sensor.



**WARNING** Pay attention to not damage the sensitive element.

8. Replace the oven heater temperature sensor in the reverse order in which it was removed.
9. Replace the oven heater baffle in the reverse order in which it was removed.
  - a. Guide the heater and temperature sensor wires through the pipes provided on the oven back wall.
  - b. Reconnect the temperature sensor and heater wires to the relevant connector **J4**, **J23**, **J24**, and **J25** on the backplane board.



- c. Plug the pipes with the insulation material.
- d. Secure the oven heater baffle to the oven wall with the four Torx screws in the tabs.

The tab adjustment screws should be loose. Carefully center the oven heater baffle inside of the oven. Be sure that the center hole in the oven heater baffle screen is centered over the blower motor shaft. The tolerance is very small. The oven heater baffle will not move left or right, nor up or down more than a couple of millimeters. However, be sure that it is centered in the oven over the blower wheel. This will ensure that whenever the GC is turned on, the blower wheel will not be touching the heater element mounted to the oven heater baffle.

10. Reinstall the back cover.
  - a. Reconnect the ground wire to the back cover terminal.
  - b. Replace the cover proceeding in the reverse order in which it was removed.
11. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
12. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
13. Close the front door of the GC. Be sure that the oven blower is running.
  - a. Verify that the heater heats, and the fan rotates without contacting the heater element.
14. Open the front door of the GC.
15. Reinstall the column.
16. Set the normal detector and GC working conditions.

## Replacing the Oven Motor

Removing the oven motor requires the removal of several parts and assemblies.



**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.


### ❖ To replace the oven motor

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
5. Remove the back cover.
  - a. Using a T20 Torxhead screwdriver remove the four screws that secure the back cover to the GC. See [Figure 208](#).

### Figure 208. Back Cover Removal

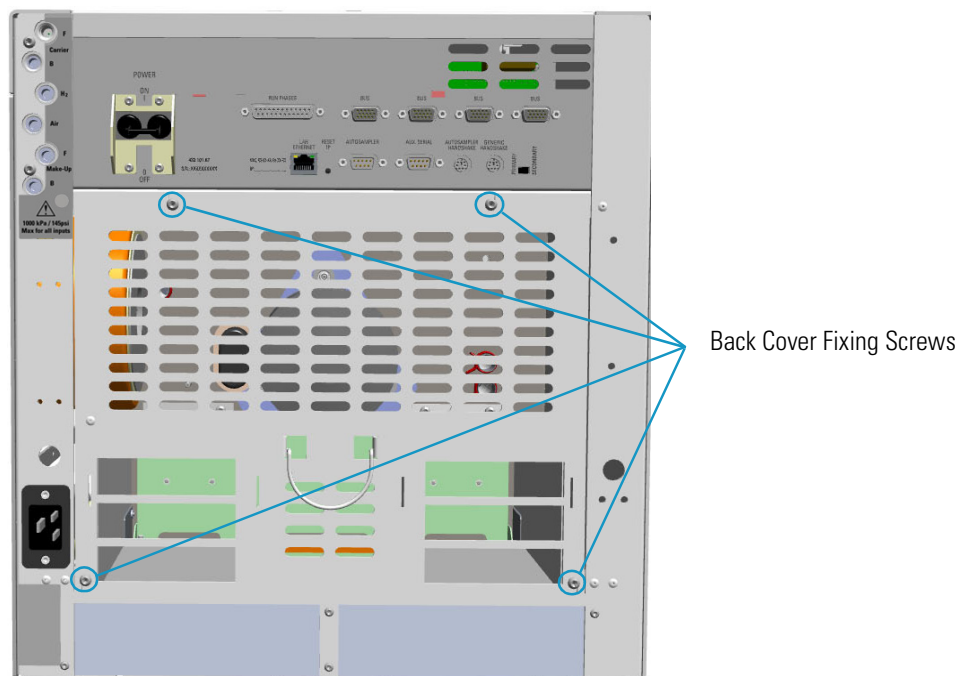

- b. Lift the cover off using the cover handle. Be aware that the ground wire is attached to the back panel.


**Note** Pay attention to the positioning of the ground wire plug, so it can be reconnected in the same way it was removed.

6. Remove the analytical column.
  - a. Open the front door of the GC.

- b. Loosen the retaining nut from the injector and detector fitting on the upper interior wall of the GC oven.
  - c. Remove the analytical column with its nut and ferrule from the bottom of the injector and the detector.
7. Remove the oven heater baffle.

The oven heater baffle is attached to the oven wall by four adjustable tabs. These tabs are attached to the oven heater baffle with four Torx screws. See [Figure 209](#).

**Figure 209.** Oven Heater Baffle

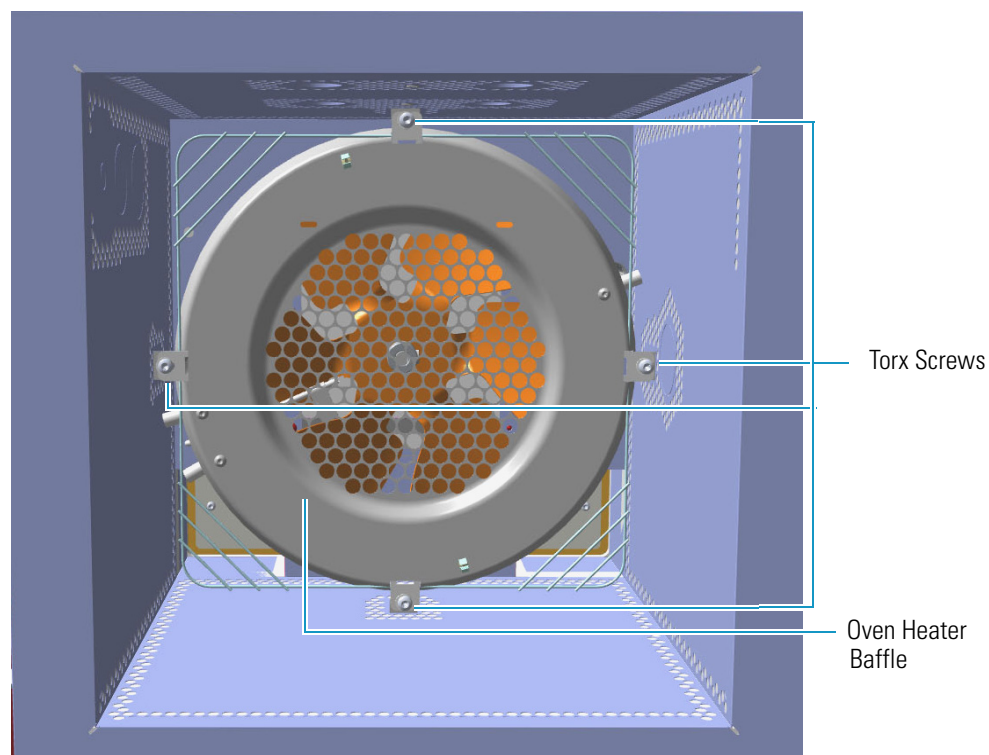


The heater and temperature sensor wires are connected to the backplane board passing through the three pipes provided on the oven back wall.

- The temperature sensor wire is plugged into the connector marked **J4** (Oven PT100)
- The heater wires are plugged into the connector **J23**, **J24**, and **J25** (Oven Heater).

The pipes are plugged by insulation material. See [Figure 210](#).

**Figure 210.** Pipes

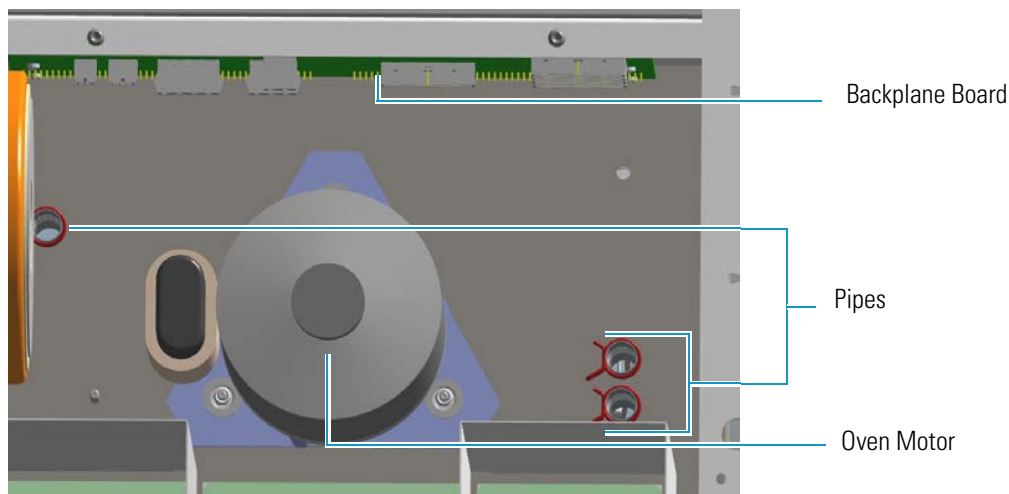

- Unplug the heater and temperature sensor wires from the relevant connector on the backplane board.
- Loosen and remove the four Torx screws that secure the tabs to the oven heater baffle. See [Figure 209](#) on [page 221](#).
- Place the heater baffle on the roof of the oven.




**IMPORTANT** During this operation, parts of the insulation material drops down into the oven. This material is reused later when a new Oven Heater Baffle is installed.

- Remove the blower fan.


The oven blower fan is attached to the shaft of the oven blower motor with an Allen screw. See [Figure 211](#).

**Figure 211.** Oven Blower Fan

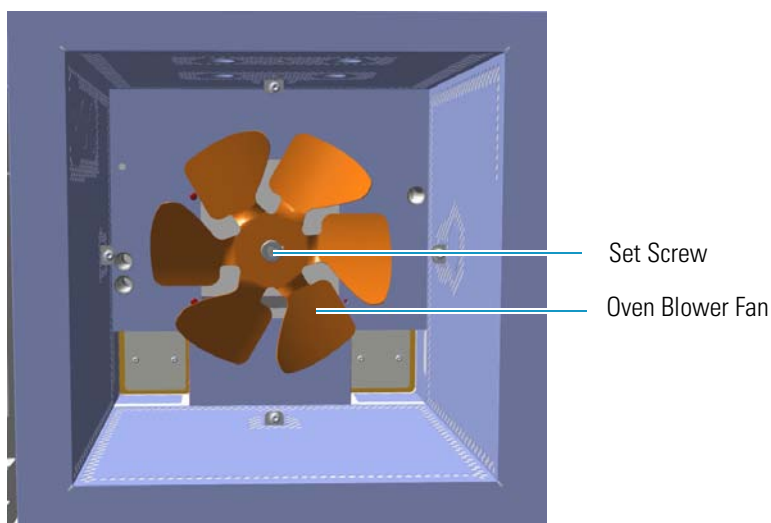

- a. Loosen the Allen screw. Carefully remove the oven fan from the blower motor shaft.
9. Remove the blower motor.
- a. The blower motor plugs into the connector marked **J7** (Motor fan) on the Backplane board. Unplug the blower motor electrical plug by pushing down on the plug tab and pulling out.
  - b. The ground strap from the motor is secured to chassis ground on a terminal just below the starting capacitor. Locate this terminal.
  - c. Remove the terminal nut using a 4-mm nut driver; or, unplug the connector. Remove the ground wire.
  - d. The blower motor is secured to the back of the oven wall with three #4 self-locking nuts. Remove these nuts using the same 4-mm nut driver. See [Figure 212](#).

**Figure 212.** Removal the Oven Blower Motor

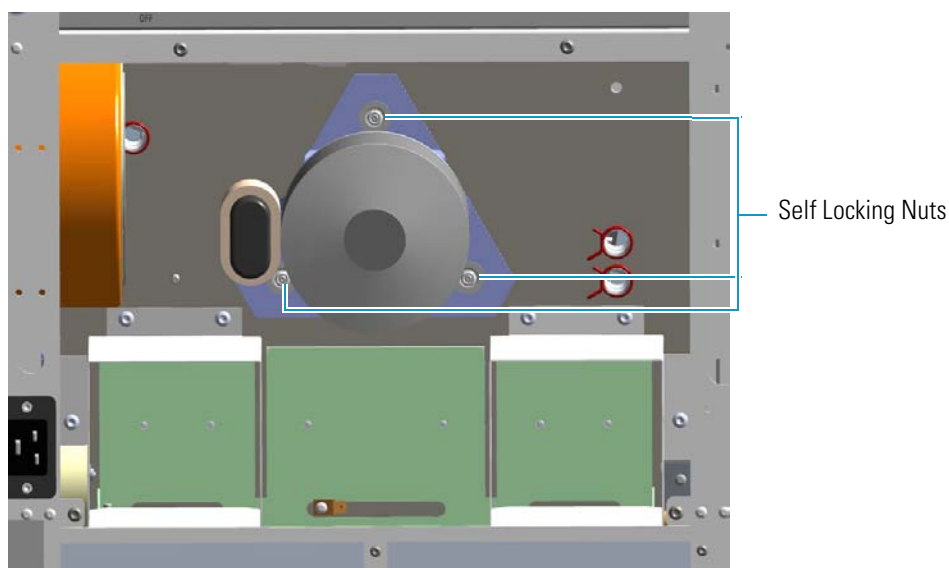

- e. Locate and note the orientation of the six shoulder washers on each side of the three grommets, before removing the blower motor.
  - f. Locate the three flat washers between the locking nut and the shoulder washers. Remove and retain these washers for reassembly.
  - g. Remove the blower motor from the oven by pulling outward.
10. Assemble the oven motor

The replacement blower motor contains the mounting brackets with the three grommets in place. The existing blower motor will have a total of six shoulder washers. Three shoulder washers are installed on each side of the grommets.

- a. Remove the all of the shoulder washers from the old blower motor.
- b. Place three of the shoulder washers over the mounting studs located on the back oven wall.

- c. Install the motor on the studs.
  - d. Place the remaining three shoulder washers over the studs, mating them with the grommets and the blower motor mounting bracket.
  - e. Place the three flat washers over the studs and against the shoulder washers.
  - f. Carefully screw the locking nuts back onto the studs, after all of the washers are in place.
  - g. Tighten the locking nuts until the grommets are flat. Loosen each locking nut by one complete turn.
  - h. Attach the ground wire to the grounding terminal on the back oven wall. Secure it with the locking nut.
  - i. Connect the blower motor electrical plug to **J7**.
11. Reassemble the remaining components on the back of the GC in the reverse order in which they were removed.
12. Reinstall the blower fan onto the shaft.
  - a. The fan should be even with the beveled edge on the motor shaft. The motor shaft should not extend more than 3 mm beyond the blower fan.
  - b. Secure the fan to the motor shaft using the M4 Allen wrench. Be sure that the setscrew is in place on the flat side of the blower motor shaft.
13. Replace the oven heater baffle in the reverse order in which it was removed.
  - a. Guide the heater and temperature sensor wires through the pipes provided on the oven back wall.
  - b. Reconnect the temperature sensor and heater wires to the relevant connector **J4**, **J23**, **J24**, and **J25** on the backplane board.
  - c. Plug the pipes with the insulation material.
  - d. Secure the oven heater baffle to the oven wall with the four Torx screws in the tabs.

The tab adjustment screws should be loose. Carefully center the oven heater baffle inside of the oven. Be sure that the center hole in the oven heater baffle screen is centered over the blower motor shaft. The tolerance is very small. The oven heater baffle will not move left or right, nor up or down more than a couple of millimeters. However, be sure that it is centered in the oven over the blower wheel. This will ensure that whenever the GC is turned on, the blower wheel will not be touching the heater element mounted to the oven heater baffle.
14. Reinstall the back cover.
  - a. Reconnect the ground wire to the back cover terminal.
  - b. Replace the cover proceeding in the reverse order in which it was removed.
15. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.

16. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
17. Close the front door of the GC. Be sure that the oven blower is running.
  - a. Verify that the heater heats, and the fan rotates without contacting the heater element.
18. Open the front door of the GC.
19. Reinstall the column.
20. Set the normal injector, detector and GC working conditions.

## Replacing the Flap Motor

To replace the oven flap motor, you must remove the complete motor flap and cooling duct assembly. The flap motor plugs into the motherboard on the GC. Therefore, you must remove the right side panel of the GC.


### ❖ To replace the flap motor

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
5. Remove the back cover.
  - a. Using a T20 Torxhead screwdriver, remove the four screws that secure the back cover to the GC. See [Figure 213](#).

**Figure 213.** Back Cover Removal

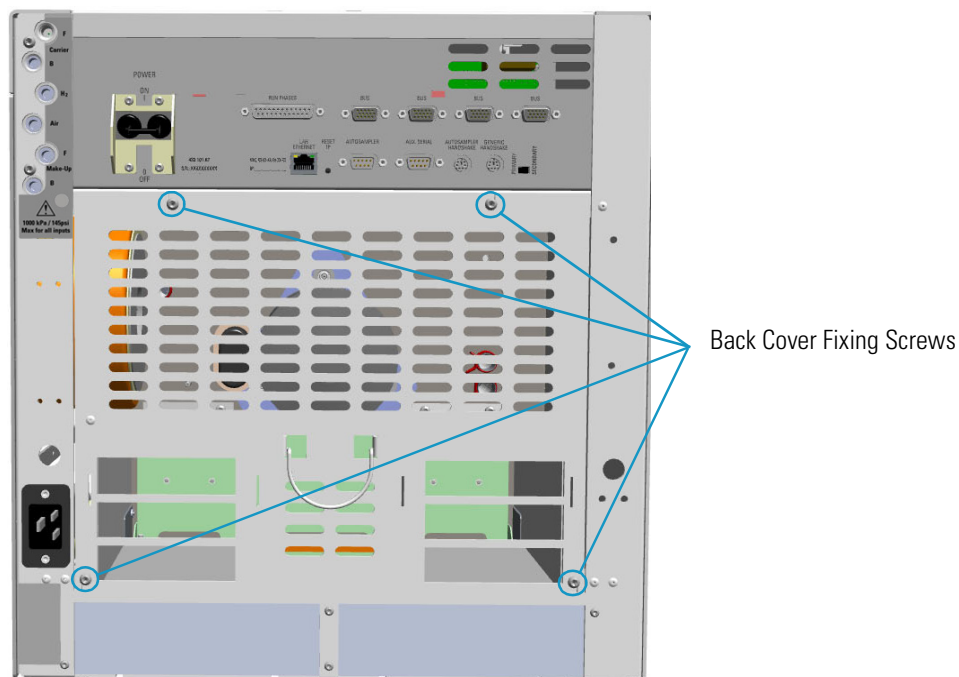


- b. Lift the cover off using the cover handle. Be aware that the ground wire is attached to the back panel.

**Note** Pay attention to the positioning of the ground wire plug, so it can be reconnected in the same way it was removed.

6. Remove the AI/AS autosampler if present.
7. Remove the right panel.
  - a. Open the front door of the GC.
  - b. Using a T20 Torxhead screwdriver remove the screw located on the upper right corner of the frontal frame.
  - c. Slide the panel towards the back of the instrument. Tilt the top of the right panel outwards.
  - d. Continue to slide the panel towards the back of the GC and lift it off.
8. Remove the cooling duct. See [Figure 214](#).




**Figure 214.** Cooling Duct Removal

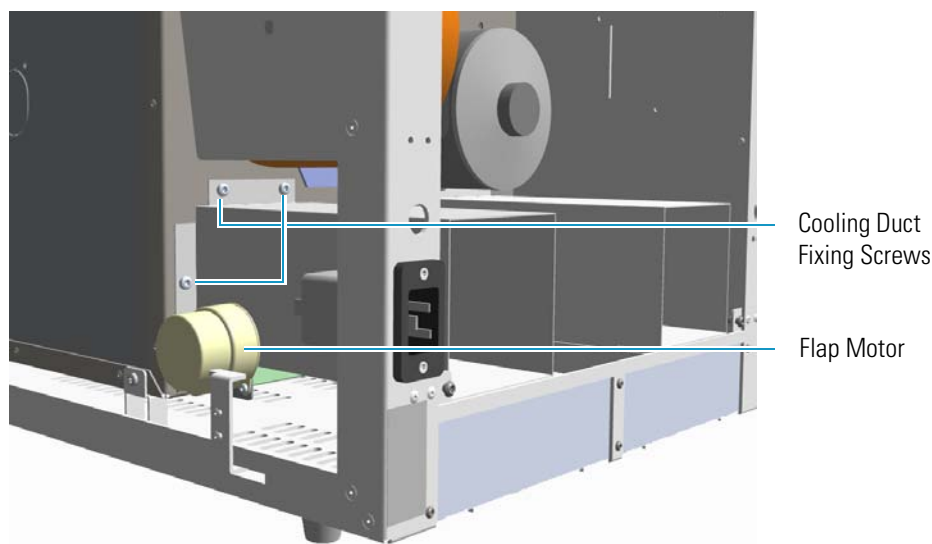

- a. Locate the three screws that secure the cooling duct to the oven wall. Remove these screws.
  - b. Remove the cooling duct from the back of the GC.
9. Remove the flap motor.
- a. Disconnect the oven flap motor electrical plug from **J10** (Flap Motor) on the backplane board.
  - b. Loosen the M4 Allen screw that secures the flap motor shaft to the railing that attaches to the flaps.
  - c. Remove the two fixing screws that secure the flap motor to its support plate. See [Figure 215](#).

**Figure 215.** Flap Motor Removal (1)

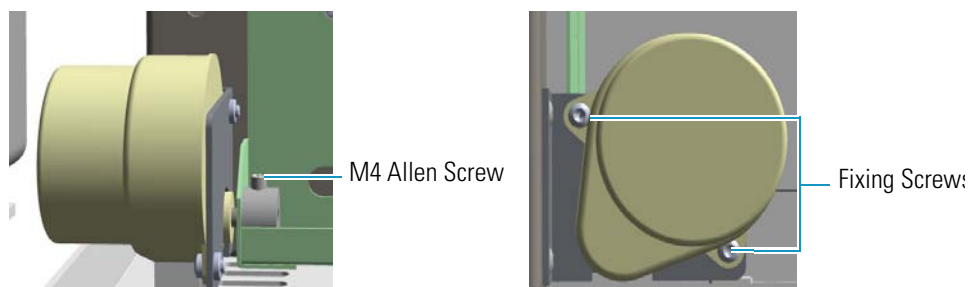

- d. Remove the flap motor paying attention to not rotate it.
10. Reinstall the new motor.
11. Align the flat side of the motor shaft with the Allen screw.
12. Reassemble the flap motor and cooling duct assembly in the reverse order in which it was removed.

13. Reconnect the power plug to connector **J10** on the backplane board.
14. Reinstall the back cover.
  - a. Reconnect the ground wire to the back cover terminal.
  - b. Replace the cover proceeding in the reverse order in which it was removed.
15. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
16. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
17. Set the normal injector, detector and GC working conditions.

## Injectors Advanced Maintenance

This chapter describes TRACE 1300/TRACE 1310 injector modules that do not require routine maintenance, but troubleshooting may indicate they need to be cleaned or replaced.

### Contents

- [Baking-out Contaminants from SSL, SSLBKF, HeS-S/SL, PTV, and PTVBKF Injectors](#)
- [Removing/Replacing an Injector Module](#)
- [Cleaning the SSL Injector Body](#)
- [Cleaning the SSLBKF Injector Body](#)
- [Cleaning the HeS-S/SL Injector Body](#)
- [Cleaning the PTV Injector Head Assembly](#)
- [Cleaning the PTVBKF Injector Head Assembly](#)

# Baking-out Contaminants from SSL, SSLBKF, HeS-S/SL, PTV, and PTVBKF Injectors

Perform this operation to bake-out contaminants from SSL, SSLBKF, HeS-S/SL, PTV, and PTVBKF injectors.


### ❖ To bake-out contaminants from the injector

1. Put the inlet into **Split** injection mode.
2. Select **Constant Flow** mode and enter the normal operating **Column flow** value.
3. Set **Split flow** to 200 mL/min.
4. Purge the column with carrier flow for at least 10 minutes before heating the oven.
5. If the column is connected to the detector, set the detector 25 °C above normal operating temperature. If the column is not attached to the detector, cap the detector fitting. In case of PDD detector, it is strongly recommended disconnecting the column from the detector.
6. Set the injector temperature to 300 °C or 25 °C above the normal operating temperature to bake out contaminants from the injector, mostly through the split vent.
7. Set the oven temperature 25 °C above the GC method final oven temperature to bake contaminants from the column. Do not exceed the column manufacturer's maximum temperature limit.
8. Bake-out for 30 minutes or until the detector baseline is free of contamination peaks.


## Removing/Replacing an Injector Module

This section provides the instructions for removing/replacing a front/back injector module.


**Figure 216.** Example of Injector Module

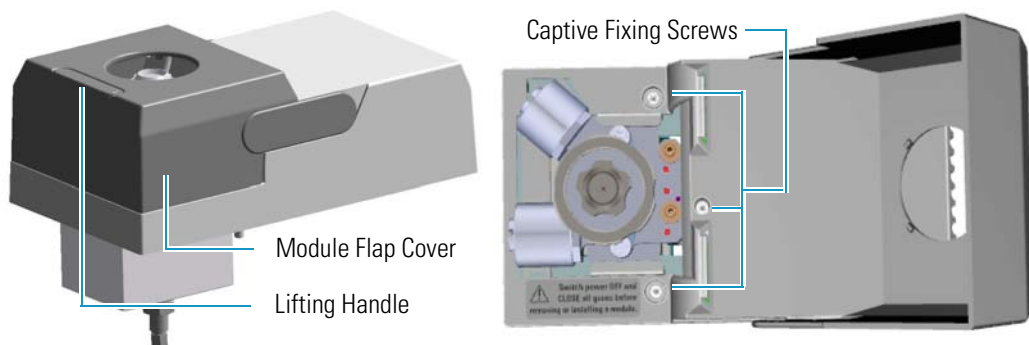

**Figure 217.** Replacing a Front Injector Module

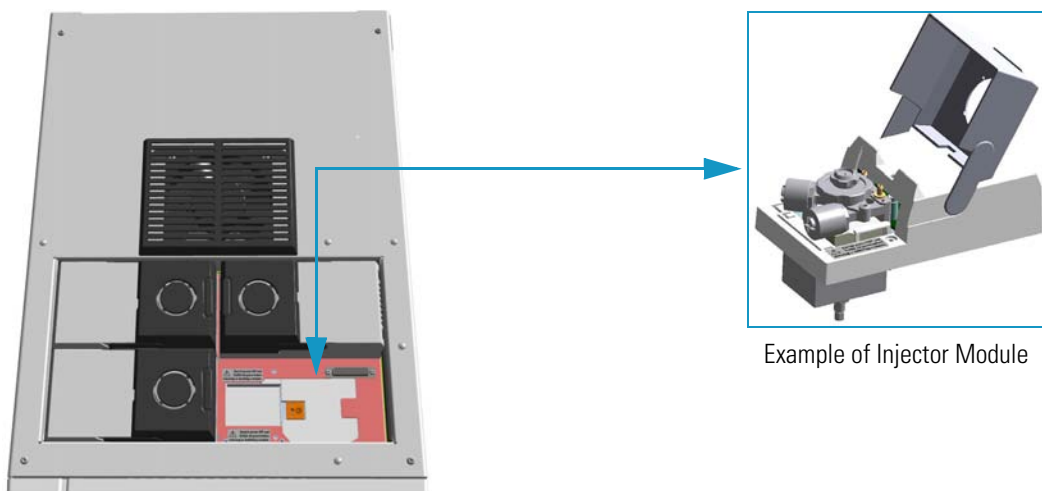

**Figure 218.** Replacing a Back Injector Module

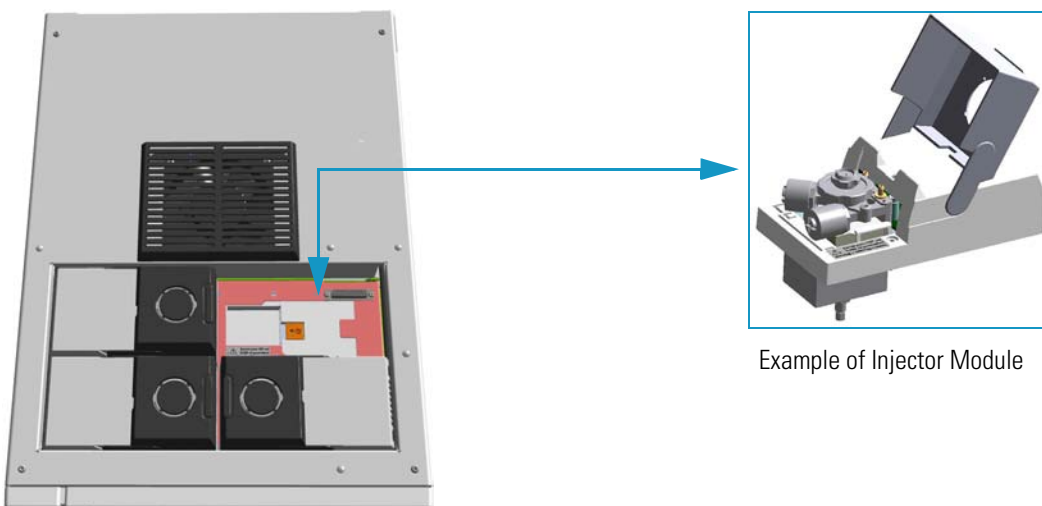


#### ❖ To remove/replace a front/back injector module

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector on the back of the GC, and from the wall outlet.
4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
5. Close the gas supplies.
6. Remove the column end from the injector.
  - a. Open the front door of the GC.
  - b. Loosen the retaining nut from the injector fitting on the upper interior wall of the GC oven.
  - c. Remove the analytical column with its nut and ferrule from the bottom of the injector.
7. Put the autosampler away if present.
  - a. If an AI/AS autosampler is present, move the sampler support to the right.
  - b. If a TriPlus or a TriPlus RSH autosampler is present, move the turret/head until enough free space is created around the module.
8. Remove the injector module.
  - a. Open the module flap cover.
  - b. Using a T20 Torxhead screwdriver, unscrew the three captive fixing screws.
  - c. Keeping the module flap cover open, lift up the module from its seat in the injector housing.
  - d. Place the module on a clean surface.

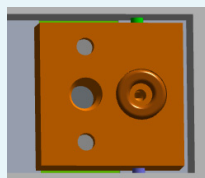

**WARNING** Make sure the o-ring is placed into its seat on the gas connection plate.

**Do not install the module if the o-ring is missing.**

9. Replace the injector module.
  - a. Open the module flap cover.
  - b. Keeping the module flap cover open, place it in its seat. Be sure to insert the 25-pin male connector, on the bottom of the module, into the 25-pin female connector on the injector seat of the injector housing.
  - c. Use a T20 Torxhead screwdriver to tighten the three captive fixing screws without overtightening.



**CAUTION** To maintain the correct alignment the screws must be tightened in turn. Tighten each screw only a small amount before moving to the next screw. Repeat until all are secure.

- d. Close the module flap cover.
  - e. Keep the plug connected to the bottom.
10. Open the gas supplies.
11. Check the gas supply for leaks.
  - a. Use a handheld electronic leak detector (Thermo Scientific GLD Pro leak detector or equivalent) to check each fitting for leaks.
  - b. If you detect a leak, tighten the connection, and retest it.
  - c. Repeat this process until all connections are leak free.
12. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
13. Power on the GC.
  - a. Plug the power cable to the AC Input connector on the back of the GC, and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
14. Pressurize the module with the carrier gas.
15. Check the module gas connections for leaks.
16. Remove the plug from the bottom.
17. Reconnect the column end to the injector, and verify the connection point.
18. Close the front door of the GC.
19. If present, move the autosampler towards the module to restore the original alignment.

## Cleaning the SSL Injector Body

Perform this operation when a more efficient cleaning of the injector body is desired or required due to contaminants.

### ❖ To clean the injector body

#### Materials needed

Ultrasonic cleaning bath

Methanol/acetone mixture 1:1

GC-grade methanol

T20 Torxhead screwdriver

Forceps or tweezers


1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the top parts of the injector. See [Figure 219](#).
  - a. Unscrew the septum cap of the injector.
  - b. Unscrew the ring nut.
  - c. Remove the septum holder/liner cap with septum from the injector body head.
7. Remove the liner. See [Figure 219](#).
  - a. Use tweezers to remove the liner with the liner seal from the injector.
8. Remove the bottom parts of the injector. See [Figure 219](#).
  - a. Unscrew the capillary column retaining nut, then remove the analytical column with its ferrule from the bottom of the injector.
  - b. Unscrew the retaining nut with the washer and the base seal.




**Figure 219.** SSL Injector Components

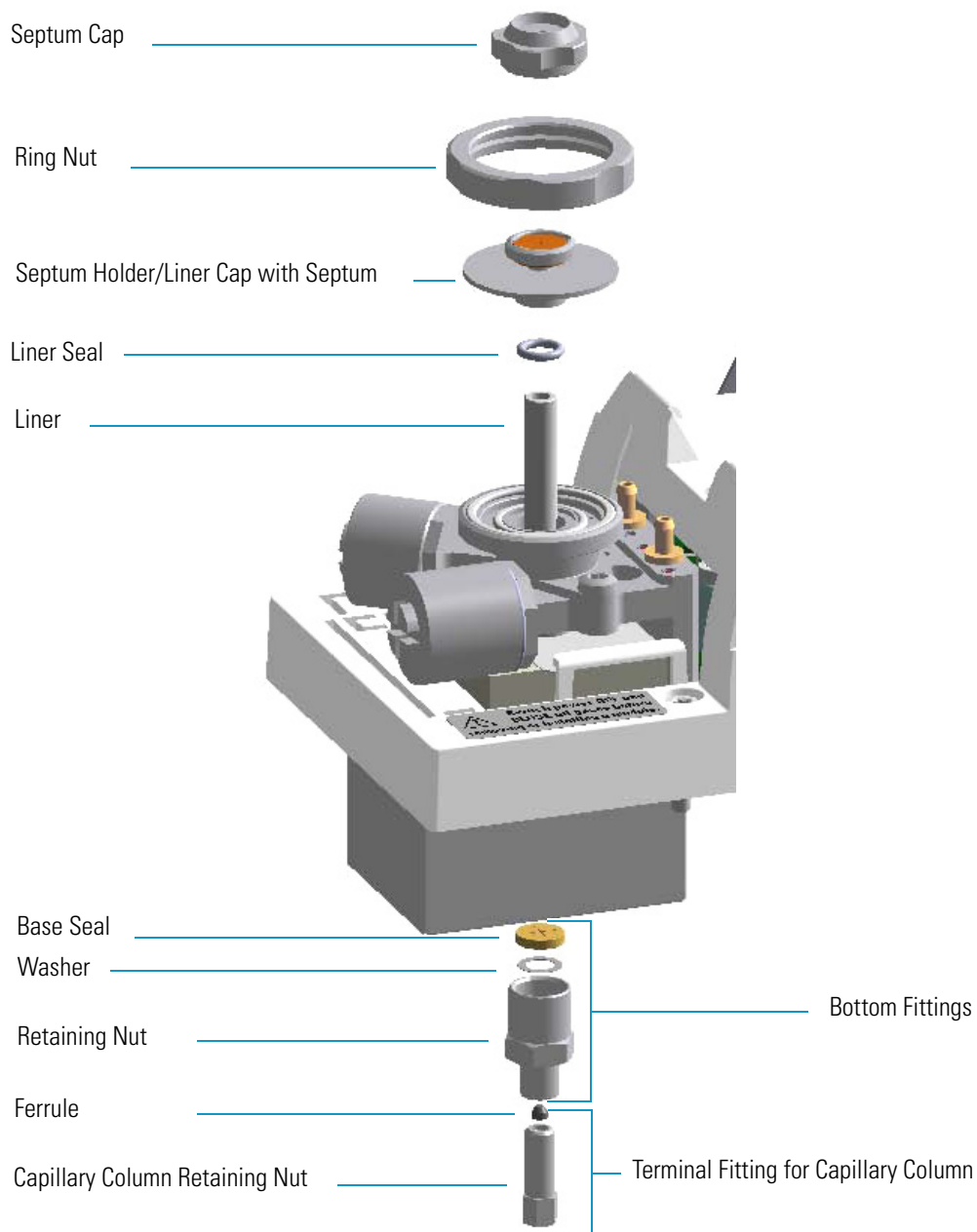

9. Remove the body head o-rings. See [Figure 220](#).
  - a. Using tweezers, remove both the internal and external body head o-rings.
  - b. Place and keep all the removed components on a clean surface.
10. Remove the injector body. See [Figure 220](#).
  - a. Using a T20 Torxhead screwdriver, undo the two injector body fixing screws, and extract the injector body from its housing.

**Note** Do not remove the carrier, split and purge lines o-rings.

**Figure 220.** SSL Injector Body Cleaning

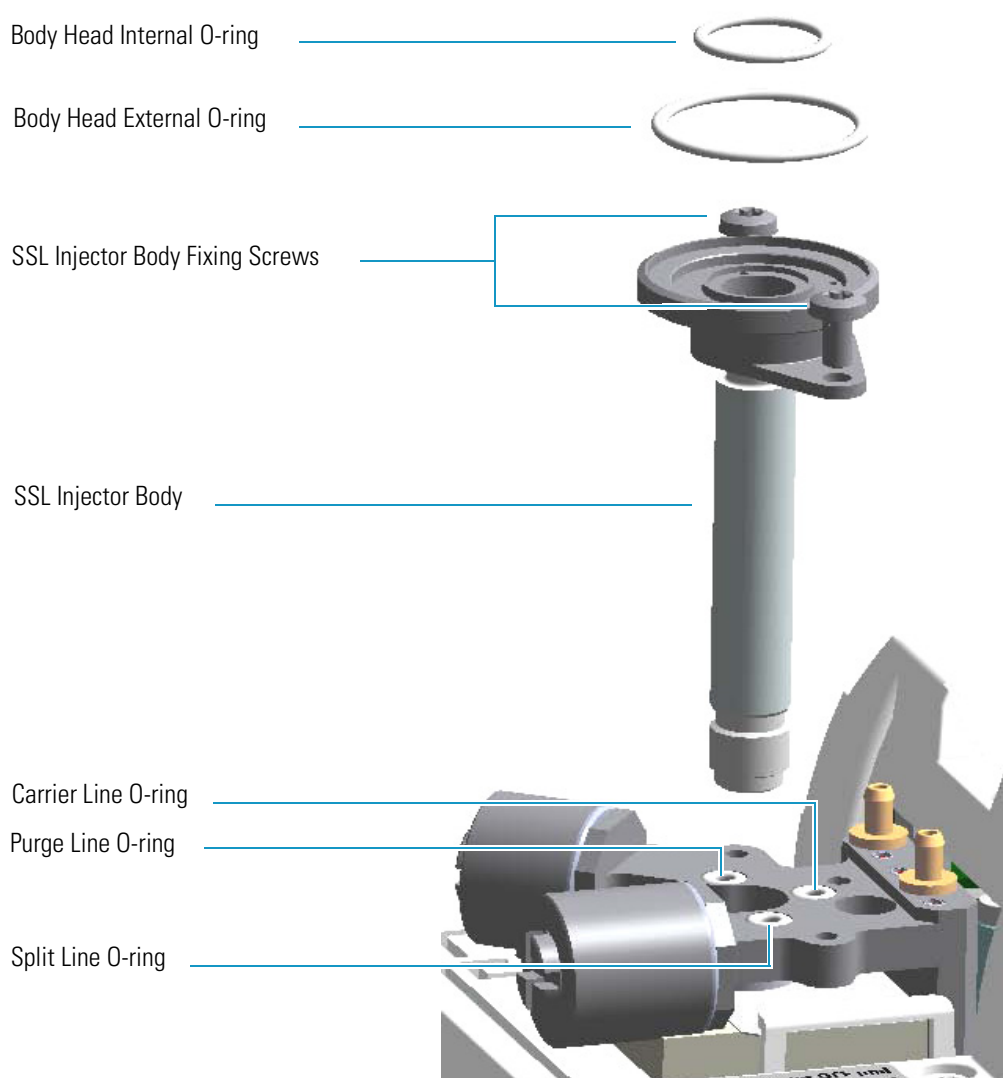

11. Clean the injector body.

- a. Put the injector body into an ultrasonic cleaner, filled with a methanol/acetone mixture (1:1).
- b. Sonicate the injector body for about half an hour.
- c. Using tweezers, remove the injector body from the bath, and dry it with an inert gas.

12. Reinstall the injector body.

- a. Reinstall and fix the injector body into its housing by screwing the two fixing screws.
- b. Using tweezers, replace both the internal and external body head o-rings.

13. Reinstall the bottom parts of the injector.
  - a. Reinstall the retaining nut with the washer and the base seal.
  - b. Reinstall the analytical column.
14. Reinstall the liner
  - a. Using tweezers, place the liner seal over the liner, insert the liner into the injector, and push it gently towards the bottom of the injector.
15. Reinstall the top parts of the injector.
  - a. Place the septum holder/liner cap with the septum on the body head of the injector. Fix them by screwing the ring nut.
  - b. Screw and tighten the injector cap to finger tight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

16. Close the module flap cover.
17. If present, move the autosampler towards the module to restore the original alignment.
18. Turn the carrier gas on.
19. Set the normal oven and injector working conditions.

## Cleaning the SSLBKF Injector Body

Perform this operation when a more efficient cleaning of the injector body is desired or required, due to contaminants.

### ❖ To clean the injector body

#### Materials needed

Ultrasonic cleaning bath

Methanol/acetone mixture 1:1

GC-grade methanol

1/8-inch -wrench

T20 Torxhead screwdriver


Forceps or tweezers

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the top parts of the injector. See [Figure 221](#).
  - a. Unscrew the septum cap of the injector.
  - b. Unscrew the ring nut.
  - c. Remove the septum holder/liner cap with septum from the injector body head.
7. Remove the liner. See [Figure 221](#).
  - a. Use tweezers to remove the liner with the liner seal from the injector.
8. Remove the bottom parts of the injector. See [Figure 221](#).
  - a. Unscrew the capillary column retaining nut, then remove the analytical column with its ferrule from the bottom of the injector.
  - b. Unscrew the retaining nut with the washer and the base seal.


**Figure 221.** SSLBKF Injector Components

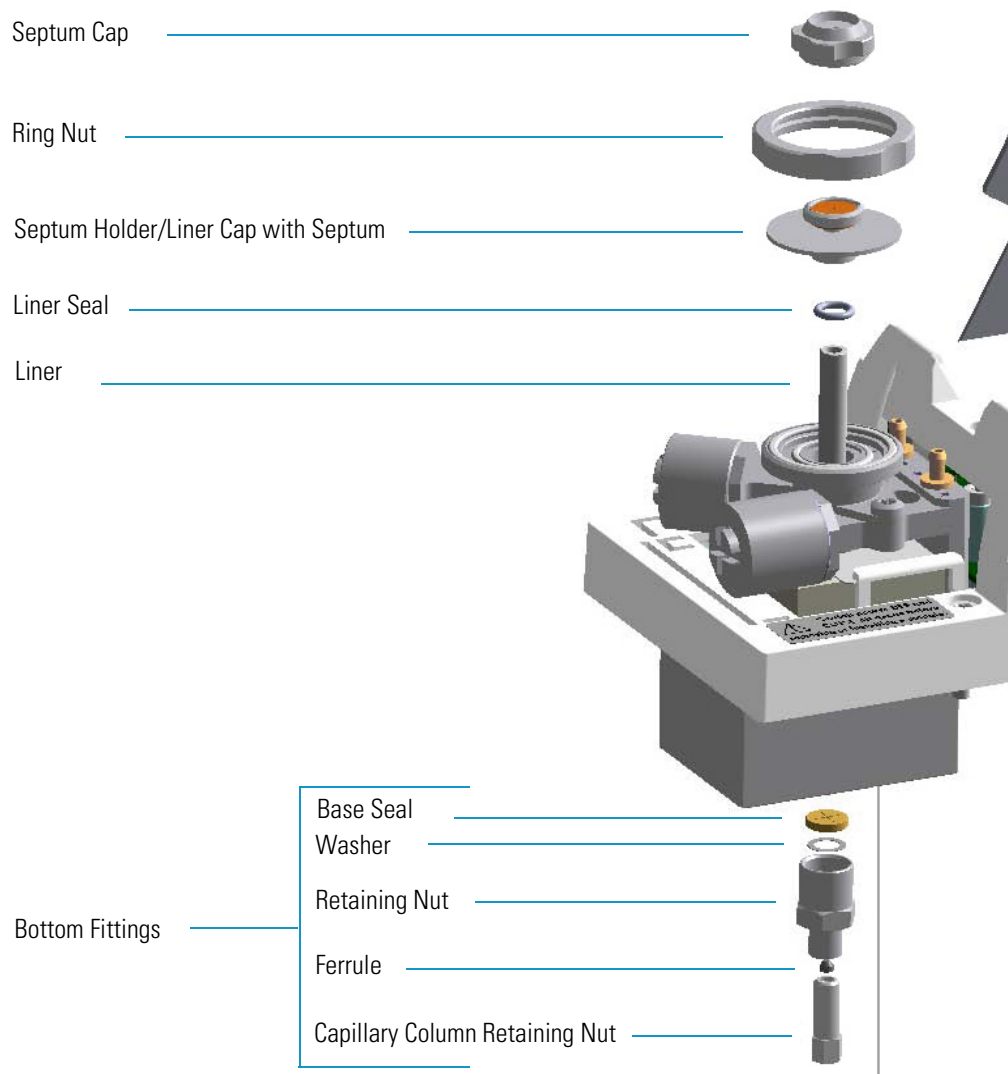

9. Remove the body head o-rings. See [Figure 222](#).
  - a. Using tweezers, remove both internal and external body head o-rings.
  - b. Place and keep all the removed components on a clean surface.
10. Remove the injector body. See [Figure 222](#).
  - a. Using a T20 Torxhead screwdriver, undo the two injector body fixing screws, and extract the injector body from its housing.

**Note** Do not remove the carrier, split, and purge lines' o-rings.

**Figure 222.** SSLBKF Injector Body Cleaning

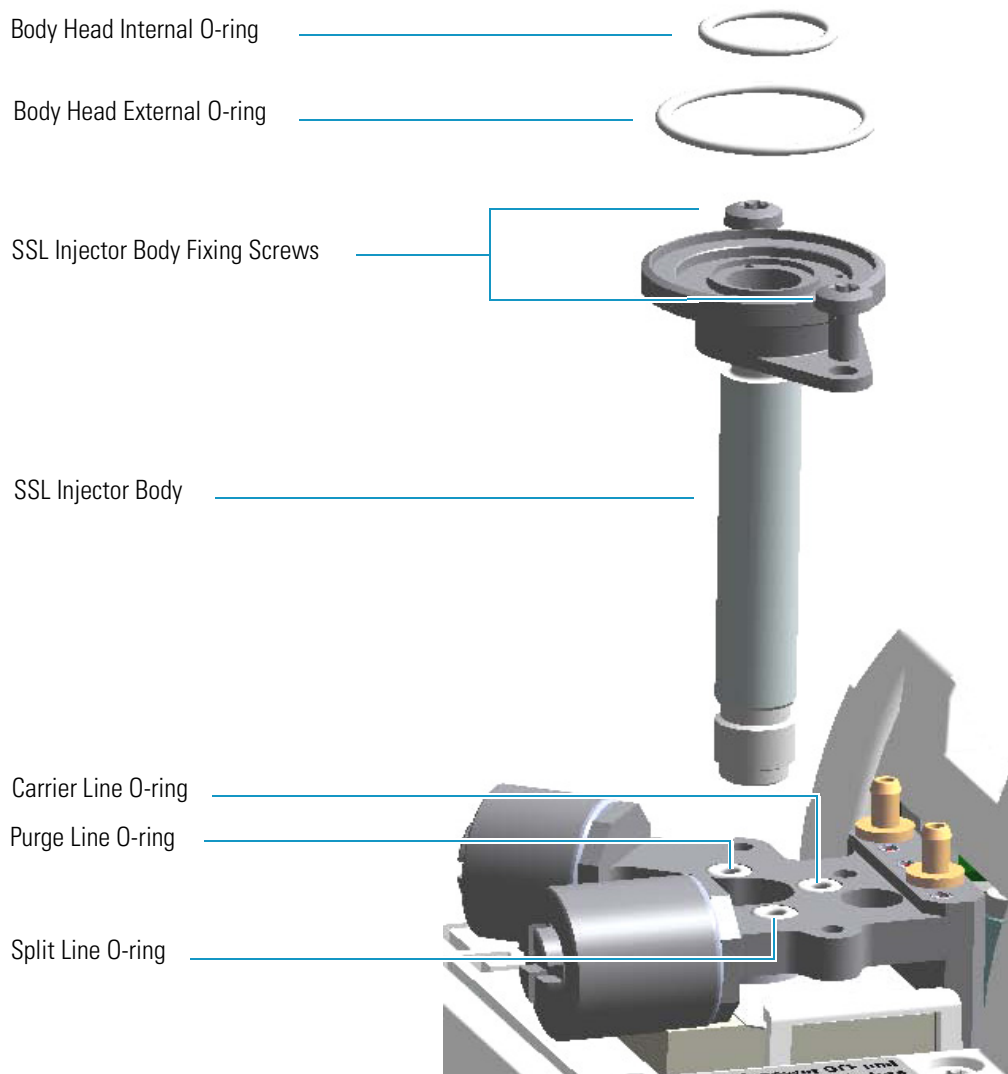

11. Clean the injector body.
  - a. Put the injector body into an ultrasonic cleaner filled with a methanol/acetone mixture (1:1).
  - b. Sonicate the injector body for about half an hour.
  - c. Using tweezers, remove the injector body from the bath, and dry it with an inert gas.
12. Reinstall the injector body.
  - a. Reinstall and fix the injector body into its housing by screwing the two fixing screws.
  - b. Using tweezers, replace both the internal and external body head o-rings.
13. Reinstall the bottom parts of the injector.
  - a. Reinstall the retaining nut with the washer and the base seal.

- b. Reinstall the analytical column.
14. Reinstall the liner
- a. Place the liner seal over the liner, then, Using tweezers, insert the liner into the injector and push it gently towards the bottom of the injector.
15. Reinstall the top parts of the injector.
- a. Place the septum holder/liner cap with the septum on the body head of the injector. Fix them by screwing the ring nut.
  - b. Screw and tighten the injector cap to fingertight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

16. Close the module flap cover.
17. If present, move the autosampler towards the module to restore the original alignment.
18. Turn the carrier gas on.
19. Set the normal oven and injector working conditions.

## Cleaning the HeS-S/SL Injector Body

Over time, contamination of the Instant Connect Helium Saver Injector Module will occur due to the deposition of cored septum particles or other material not captured by the glass wool of the injection port liner. In this case, the injector insert should be removed and cleaned according to the following procedure.

### ❖ To clean the Instant Connect Helium Saver Injector Module injector body

#### Materials needed

Ultrasonic cleaning bath

Methanol/acetone mixture 1:1

GC-grade methanol

Methylene chloride

Hexane

T20 Torxhead screwdriver

Forceps or tweezers


1. Put the GC in standby condition.
2. Cool the GC oven, injector, transfer line, and ion source.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Vent the mass spectrometer and set the inlet flow rate (nitrogen) to **Off**. Keep the helium enabled and pressurized as usual.
4. Put the autosampler away if present.
5. Open the module flap covers.
6. Remove the top parts of the injector. See [Figure 223](#).
  - a. Unscrew the septum cap of the injector.
  - b. Unscrew the ring nut.
  - c. Remove the septum holder/liner cap with septum from the injector body head.
7. Remove the liner. See [Figure 223](#).
  - a. Use tweezers to remove the liner with the liner seal from the injector.
8. Remove the analytical column.
  - a. Unscrew the capillary column retaining nut, then remove the analytical column with its ferrule from the bottom of the injector.




**Figure 223.** Instant Connect Helium Saver Injector Module: Cleaning the Injector Body (1)

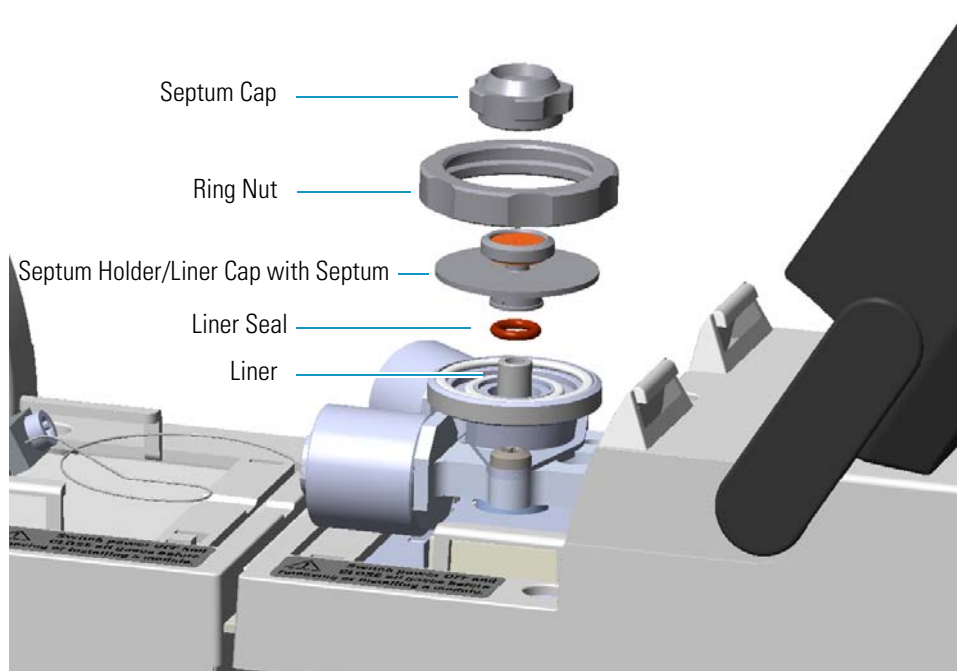

9. Disconnect the helium transfer tube See [Figure 224](#).

- a. Loosen the captive screw of the helium transfer tube and remove the tube from the gas delivery block.

**Figure 224.** Instant Connect Helium Saver Injector Module: Cleaning the Injector Body (2)

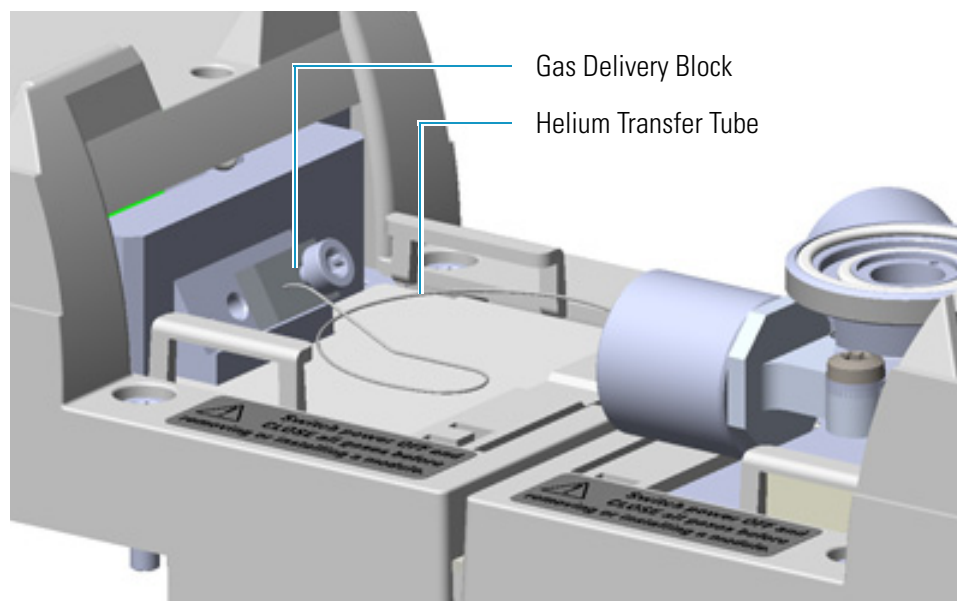


10. Remove the body head o-rings. See [Figure 225](#).

- a. Using tweezers, remove both the internal and external body head o-rings.

- b. Place and keep all the removed components on a clean surface.
11. Remove the injector body. See [Figure 225](#).
  - a. Using a T20 Torxhead screwdriver, undo the two injector body fixing screws, and extract the injector body from its housing.

**Note** Do not remove the carrier, split and purge lines o-rings.

**Figure 225.** Instant Connect Helium Saver Injector Module: Cleaning the Injector Body (3)

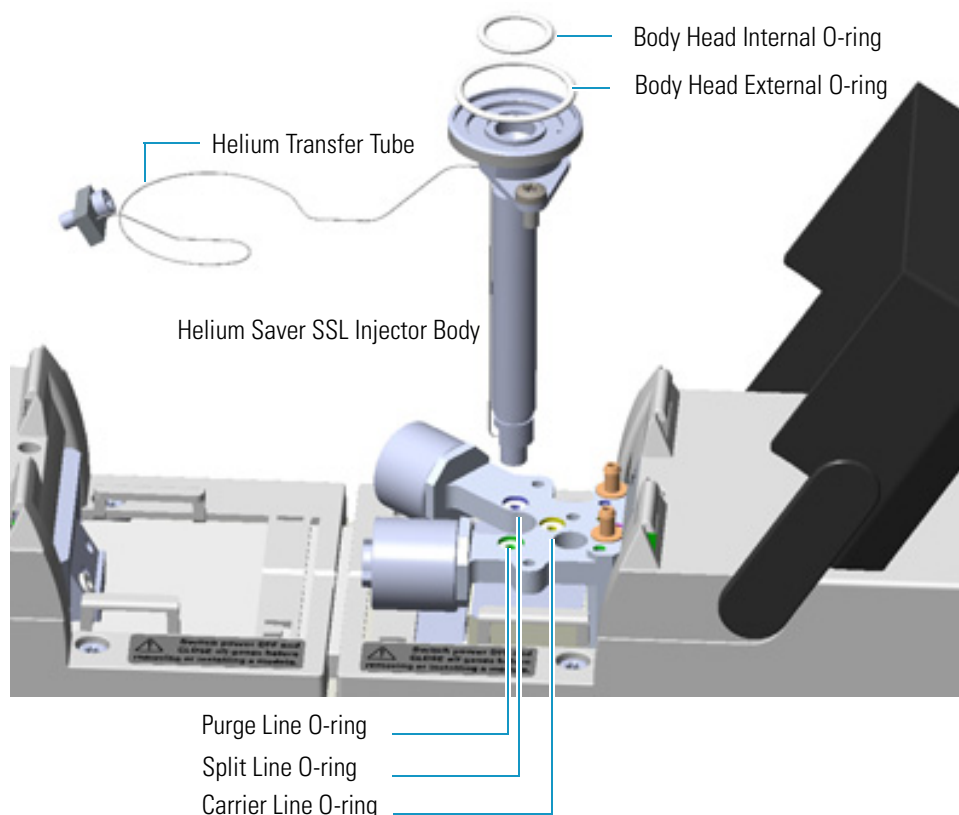

12. Clean the injector body.
  - a. Ultrasonically clean the injector insert using a warm 1% Liquinox™ solution (or equivalent).
  - b. Thoroughly rinse the insert then ultrasonically clean in chromatographic grade methanol or acetone followed by solvents of lower polarity such as methylene chloride and hexane. Limit ultra-sonication to 5 min in each solvent.
  - c. Blow dry the insert using high purity gas (do not use compressed house air as it contains residual oils from the compressor) then assemble in the reverse order of removal.

**Note** Do not use abrasives, cleaning wires, or brushes on the inlet insert as these will damage the passivation treatment layers. Be especially careful not to bend the short segment of tubing at the base of the insert interior. Nothing should be inserted into the bore of the insert other than glass liners.

13. Reinstall the injector body.
  - a. Reinstall and fix the injector body into its housing by screwing the two fixing screws.
  - b. Using tweezers, replace both the internal and external body head o-rings.
14. Reconnect the additional helium carrier gas line mating block.
15. Reinstall the liner.
  - a. Using tweezers, place the liner seal over the liner, insert the liner into the injector, and push it gently towards the bottom of the injector.
16. Reinstall the analytical column. Refer to the *TRACE 1300/1310 GC Hardware Manual*.
17. Reinstall the top parts of the injector.
  - a. Place the septum holder/liner cap with the septum on the body head of the injector. Fix them by screwing the ring nut.
  - b. Screw and tighten the injector cap to finger tight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

18. Close the module flap covers.
19. If present, move the autosampler towards the module to restore the original alignment.
20. Wait two or three minutes for helium to purge the inlet.
21. Turn the inlet flow back on and enable the inlet heater.
22. Set the normal injector, detector, and GC working conditions.

## Cleaning the PTV Injector Head Assembly

Perform this operation when cleaning the injector head assembly is desired or required, due to contaminants.



**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.

### ❖ To clean the injector head assembly

#### Materials needed

|                                       |
|---------------------------------------|
| Ultrasonic cleaning bath              |
| GC-grade Methanol/acetone mixture 1:1 |
| 1/8-inch -wrench                      |
| 5.5 x 25 mm Slotted Stubby Driver     |
| T10 Torxhead screwdriver              |
| Forceps or tweezers                   |

1. Put the GC in standby condition.
2. Cool the oven, injectors, and detectors to room temperature.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Power off the GC.
  - a. Pull down on the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector on the back of the GC, and from the wall outlet.
4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
5. Close the gas supplies.
6. Put the autosampler away if present.
7. Remove the top parts of the injector. See [Figure 226](#).
  - a. Unscrew the septum cap of the injector.
  - b. Using tweezers, remove the septum from the top of the injector head assembly.
8. Remove the liner. See [Figure 226](#).
  - a. Using the slotted stubby driver provided, unscrew and remove the liner cap.

- b. Use tweezers to remove the liner with the liner seal (o-ring) from the injector.

**CAUTION** Be careful not to break the glass liner when removing it. Glass splinters might fall into the lower part of the vaporization chamber. If the glass liner breaks, follow the instructions for [“Replacing the PTV Broken Liner”](#) on [page 123](#).

**Figure 226.** PTV Top Parts and Liner Removal

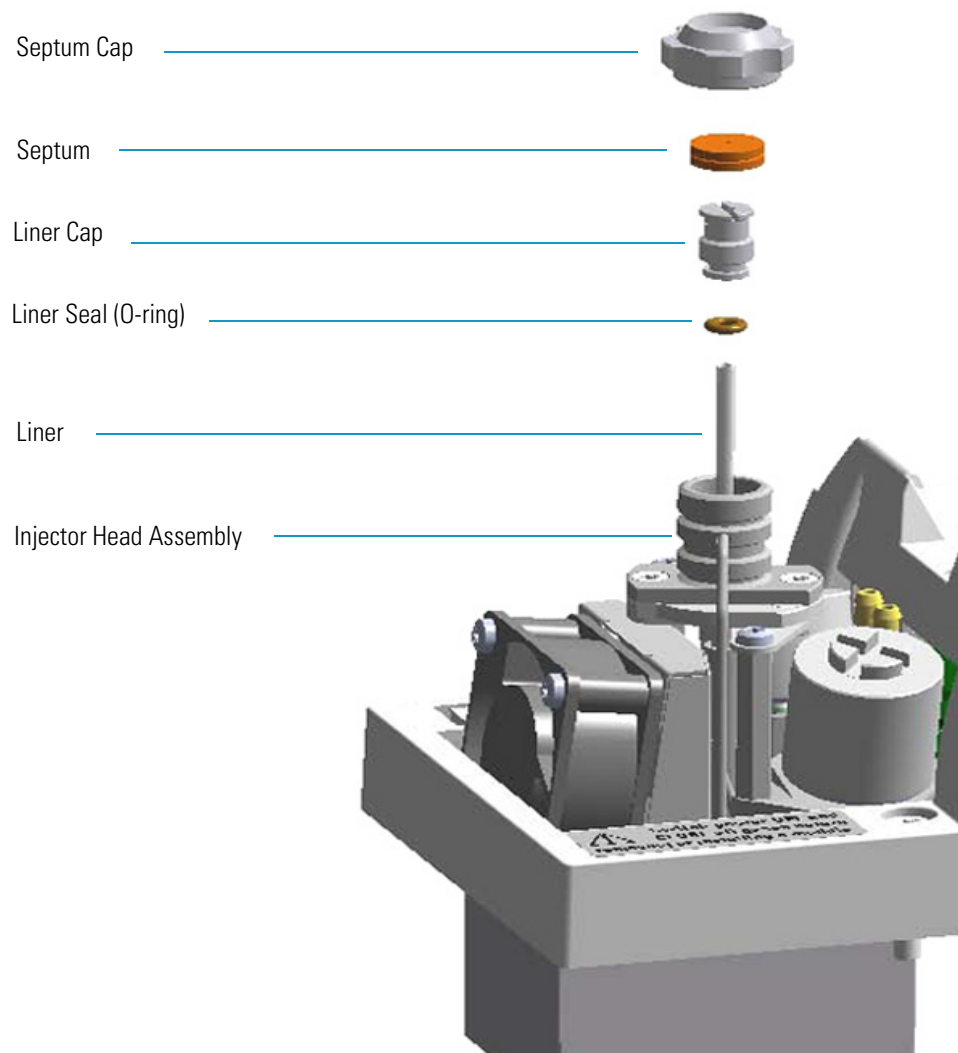


9. Remove the fan assembly.

The fan is fixed on a support. Do not remove the fan from its support but remove the assembly proceeding as follows:


- a. Using tweezers and a T10 Torxhead screwdriver, unscrew the left fixing screw of the support, then unscrew the right fixing screw. See [Figure 227](#).

**Figure 227.** PTV Module: Fan Assembly

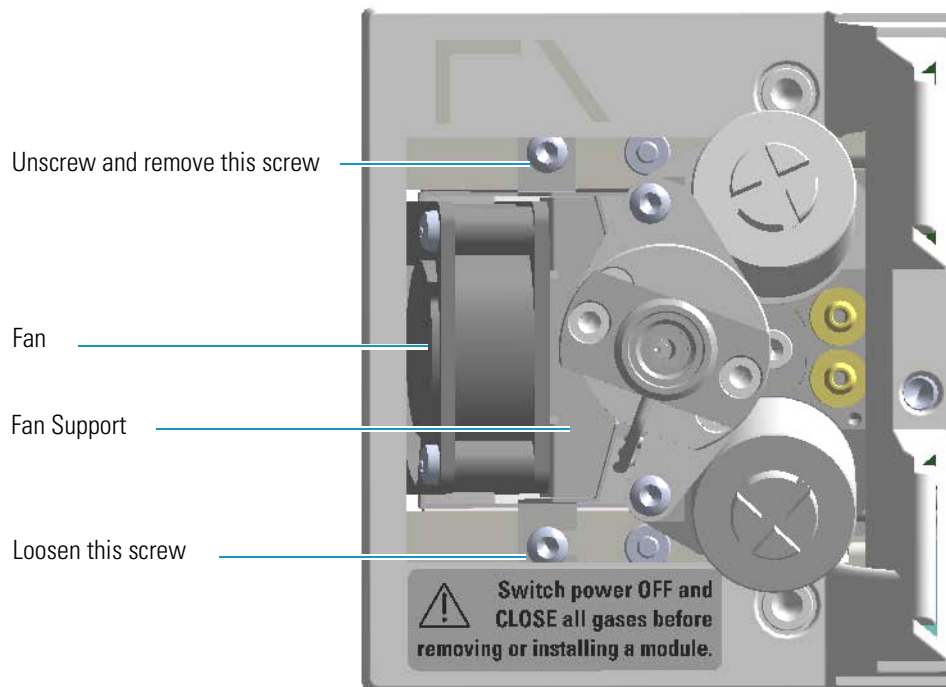


- b. Carefully turn the fan assembly slightly to access the slot on the right base of the support from the right fixing screw. See [Figure 228](#).

**Figure 228.** PTV Module: Fan Assembly Removal (1)

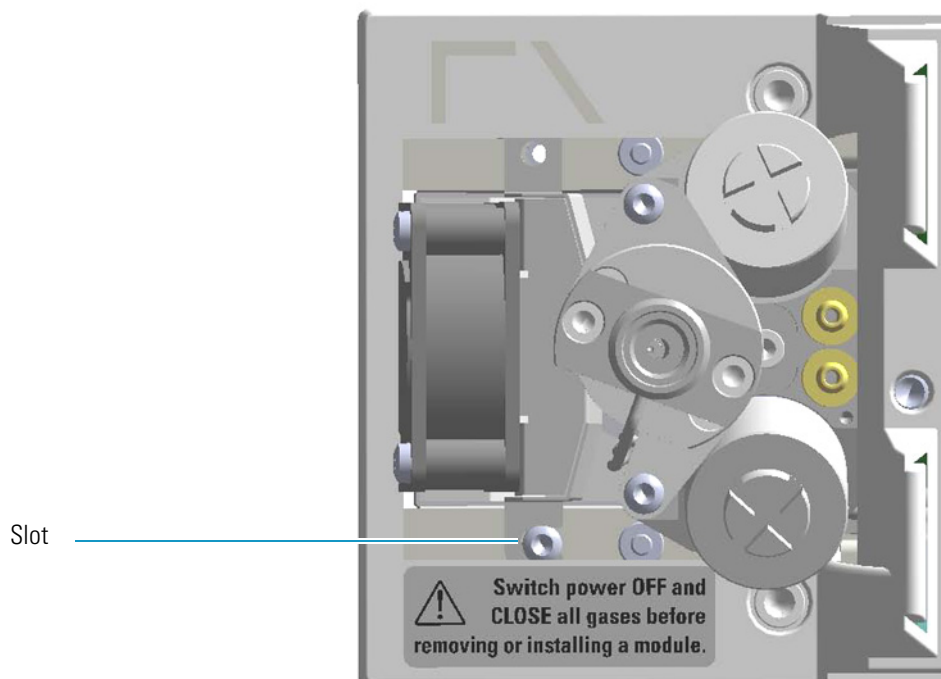

- c. Carefully extract the fan assembly from its seat, and place it sideways to the module. Be sure to not damage the two wires. See [Figure 229](#).

**Figure 229.** PTV Module: Fan Assembly Removal (2)

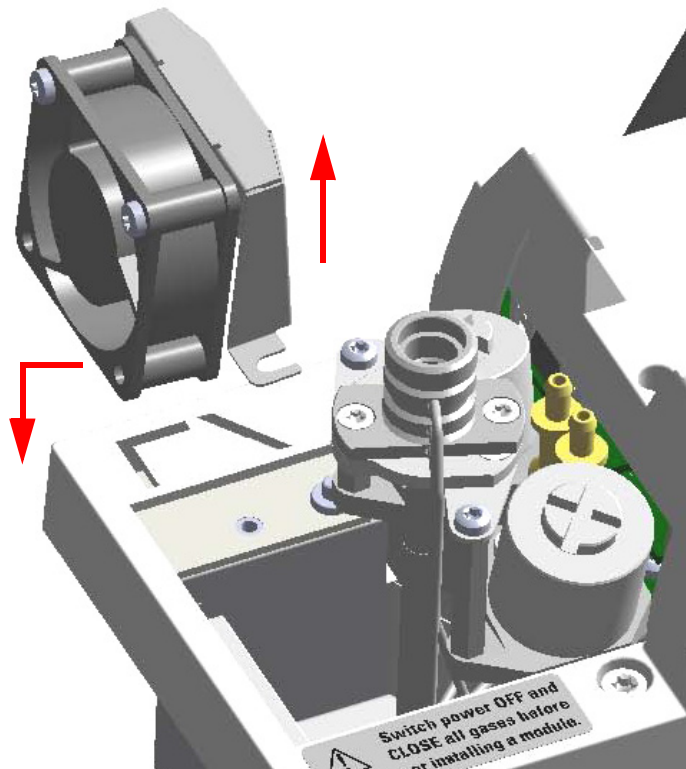


10. Remove the injector head.

- a. Using the T10 Torxhead screwdriver, unscrew the two fixing screws, and the captive screw that respectively fix the top and the bottom of the injector head assembly. See [Figure 230](#).

**Figure 230.** PTV Injector Head Assembly


Top Injector  
Head Assembly  
Fixing Screws

Bottom Injector  
Head Assembly  
Captive Screw

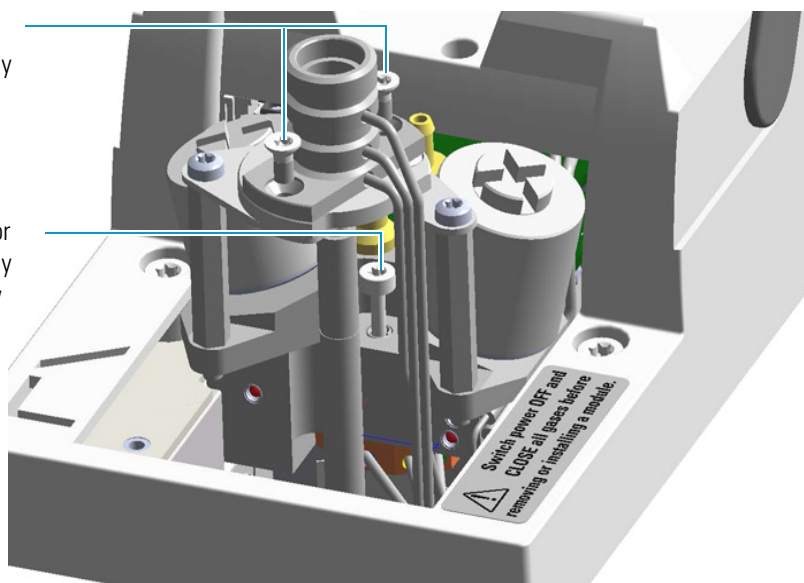

- b. Holding the top of the injector head assembly with your thumb, index, and middle fingers, turn the assembly clockwise. See [Figure 231](#).


**Figure 231.** PTV Injector Head Assembly Removal (1)

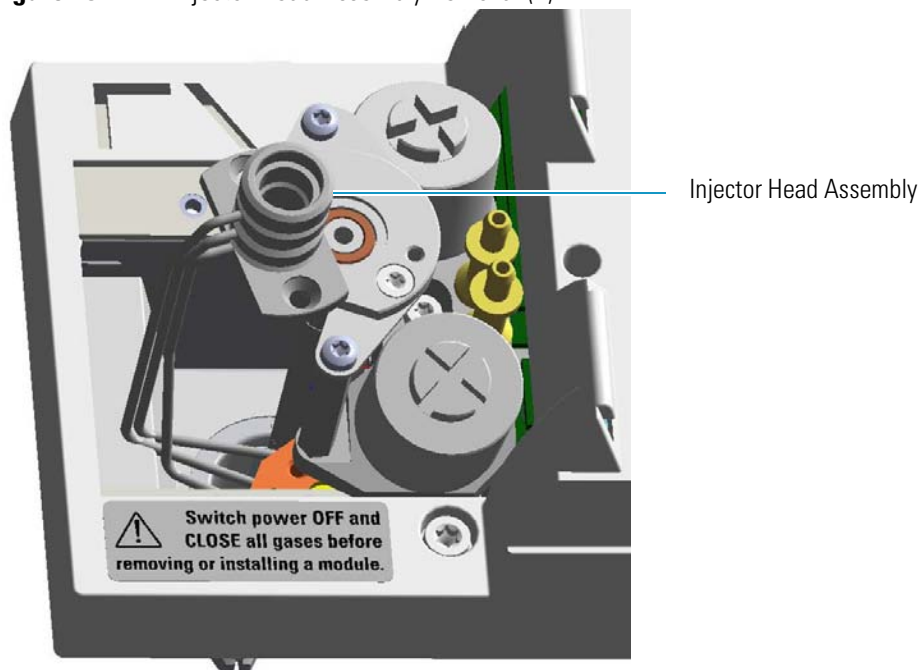

- c. Extract the assembly sideways from its housing. See [Figure 232](#).

**Note** Do not remove the carrier, split, and purge lines' o-rings on the plate.

**Figure 232.** PTV Injector Head Assembly Removal (2)

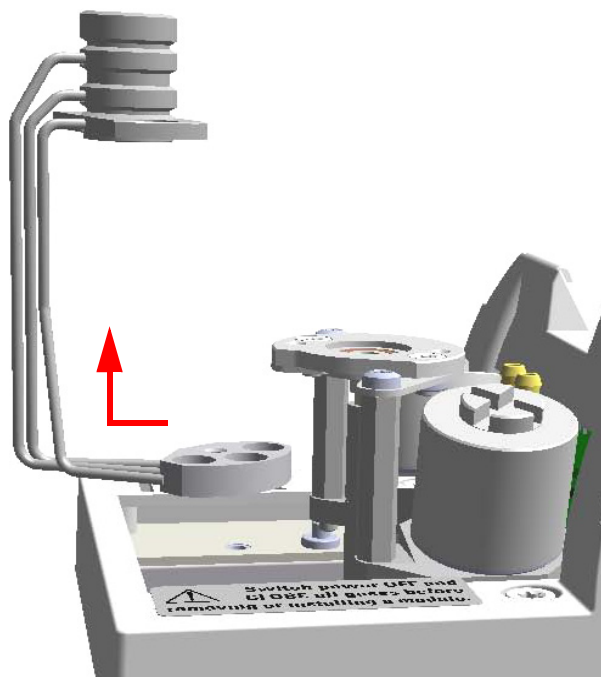



11. Clean the injector head.
  - a. Put the injector head assembly into an ultrasonic cleaner filled with a methanol/acetone mixture (1:1).
  - b. Sonicate the injector head assembly for about half an hour.
  - c. Using tweezers, remove the injector head assembly from the bath, and dry it with compressed clean air.
12. Reinstall the injector head.
  - a. Carefully reinstall the injector head assembly into its housing proceeding in the reverse order in which it was removed.
  - b. Screw the two fixing screws and the captive screw fixing the top and the bottom of the injector head assembly respectively.
13. Remount the fan assembly.
  - a. Carefully reinsert the fan assembly in its seat into the module.
  - b. Move the assembly for inserting the slot on the right base of the support slightly under the head of the right fixing screw, and for aligning the hole on the left base with the corresponding fixing hole.
  - c. Using tweezers, reinsert the left fixing screw previously removed.
  - d. Using a T10 Torxhead screwdriver tighten both the left and right fixing screws.
14. Reinsert the liner.
  - a. Place the liner seal over the liner.
  - b. Using tweezers, insert the liner into the injector, and push it gently towards the bottom of the injector.
  - c. Reinsert and screw the liner cap using the slotted stubby driver.
15. Reinstall the upper parts of the injector.
  - a. Place the septum into the injector head assembly.



**CAUTION** We suggest replacing the septum with a new one.

- b. Screw and tighten the septum cap to fingertight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

16. Close the module flap cover.
17. If present, move the autosampler towards the module to restore the original alignment.

18. Open the gas supplies.
19. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
20. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
21. Set the normal oven, injectors, and detectors working conditions.

## Cleaning the PTVBKF Injector Head Assembly

Perform this operation when cleaning the injector head assembly is desired or required, due to contaminants.



**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Field Service Engineers.

### ❖ To clean the injector head assembly

#### Materials needed

|                                       |
|---------------------------------------|
| Ultrasonic cleaning bath              |
| GC-grade Methanol/acetone mixture 1:1 |
| 1/8-inch -wrench                      |
| 5.5 x 25 mm Slotted Stubby Driver     |
| T10 Torxhead screwdriver              |
| Forceps or tweezers                   |

1. Put the GC in standby condition.
2. Cool the oven, injectors, and detectors to room temperature.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
5. Close the gas supplies.
6. Put the autosampler away if present.
7. Open the module flap cover.
8. Remove the top parts of the injector. See [Figure 233](#).
  - a. Unscrew the septum cap of the injector.
  - b. Using tweezers, remove the septum from the top of the injector head assembly.
9. Remove the liner. See [Figure 233](#).

- a. Using the slotted stubby driver provided, unscrew and remove the liner cap.
- b. Use tweezers to remove the liner with the liner seal (o-ring) from the injector.

**CAUTION** Be careful not to break the glass liner when removing it. Glass splinters might fall into the lower part of the vaporization chamber. If the glass liner breaks, follow the instructions for [“Replacing the PTVBKF Broken Liner”](#) on [page 133](#).

**Figure 233.** PTVBKF Top Parts and Liner Removal

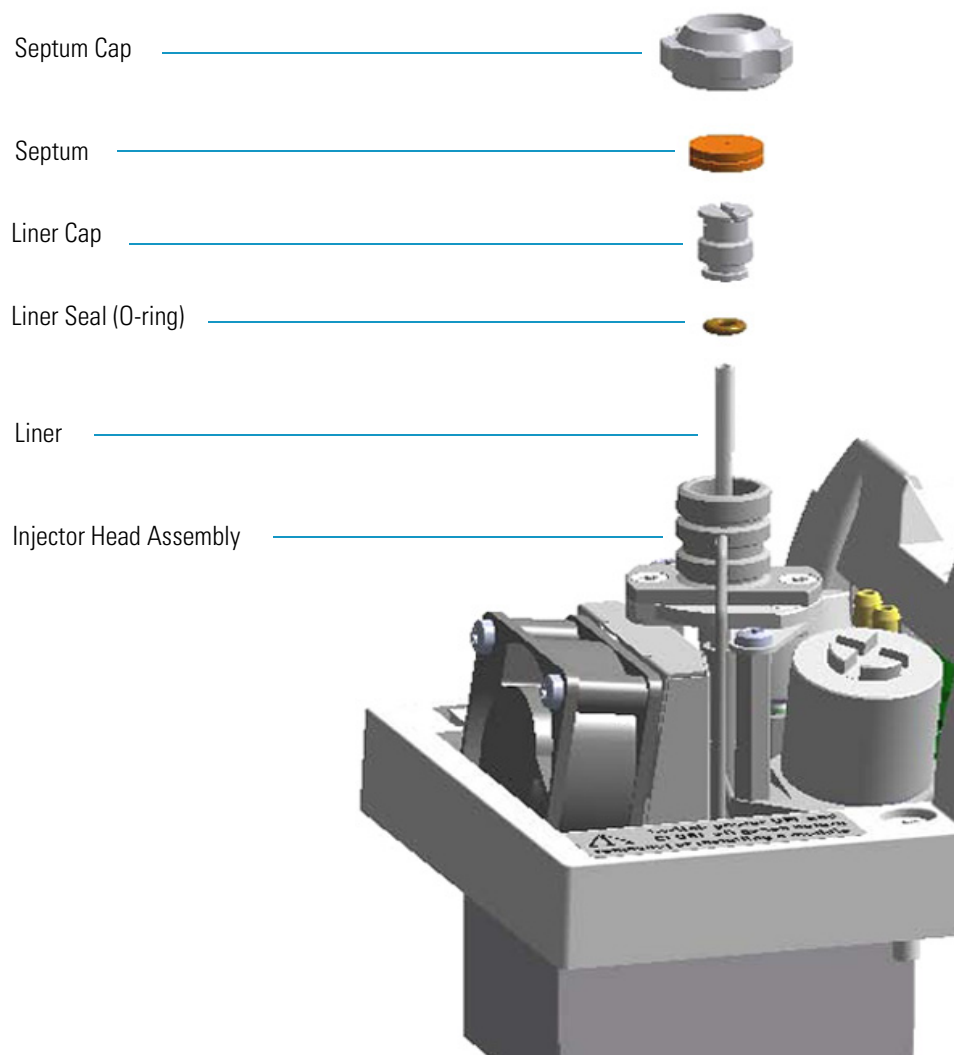


10. Remove the fan assembly.

The fan is fixed on a support. Do not remove the fan from its support, but remove the assembly by proceeding as follows:


- a. Using tweezers, and a T10 Torxhead screwdriver, unscrew and remove the left fixing screw of the support, then unscrew and remove the right fixing screw. See [Figure 234](#).

**Figure 234.** PTVBKF Module: Fan Assembly

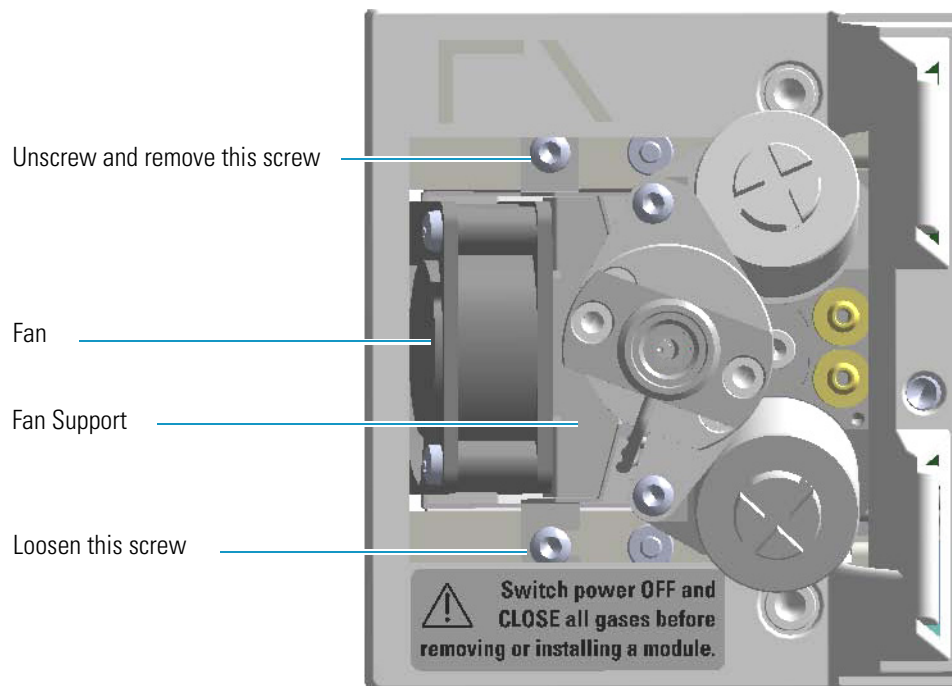


- b. Carefully turn the fan assembly slightly to access the slot on the right base of the support from the right fixing screw. See [Figure 235](#).

**Figure 235.** PTVBKF Module: Fan Assembly Removal (1)

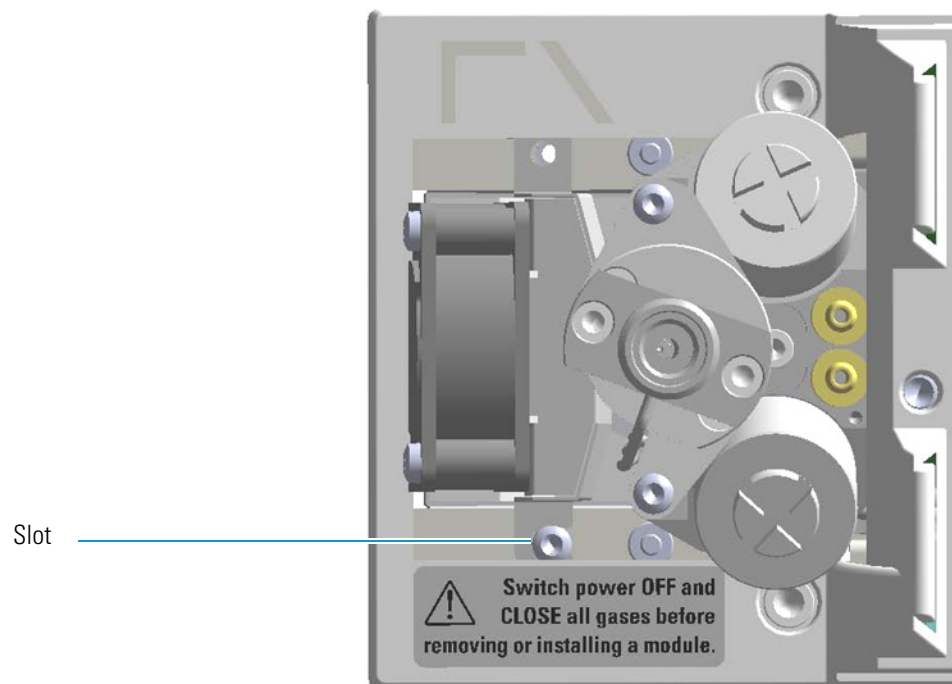

- c. Carefully extract the fan assembly from its seat, and place it sideways to the module. Be sure to not damage the two wires. See [Figure 236](#).

**Figure 236.** PTVBKF Module: Fan Assembly Removal (2)

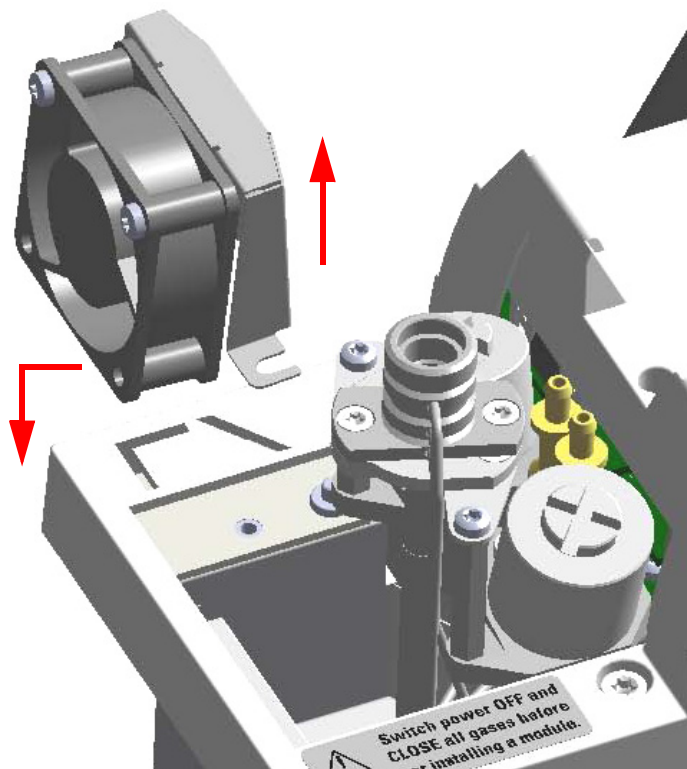


11. Remove and clean the injector head.

- a. Using the T10 Torxhead screwdriver, unscrew the two fixing screws, and the captive screw that respectively fix the top and the bottom of the injector head assembly. See [Figure 237](#).

**Figure 237.** PTVBKF Injector Head Assembly


Top Injector  
Head Assembly  
Fixing Screws

Bottom Injector  
Head Assembly  
Captive Screw

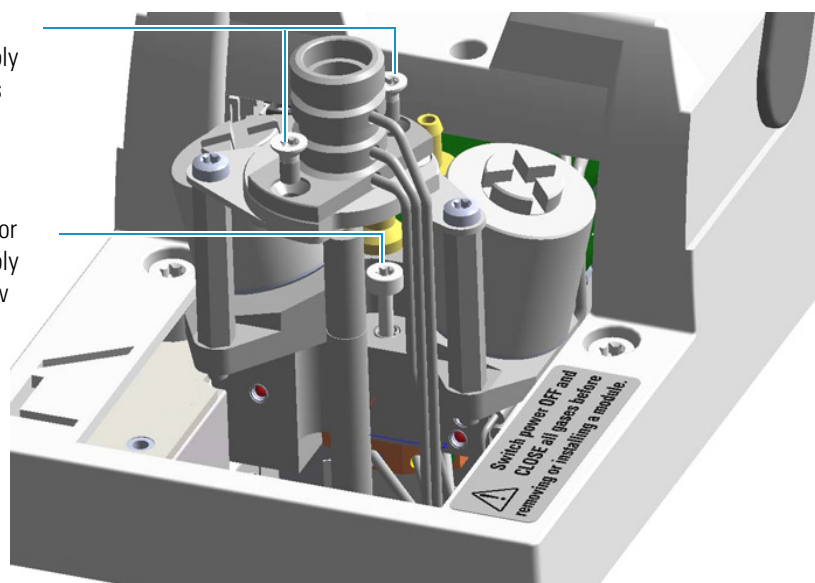

- b. Holding the top of the injector head assembly with your thumb, index, and middle fingers, turn the assembly clockwise. See [Figure 238](#).


**Figure 238.** PTV Injector Head Assembly Removal (1)

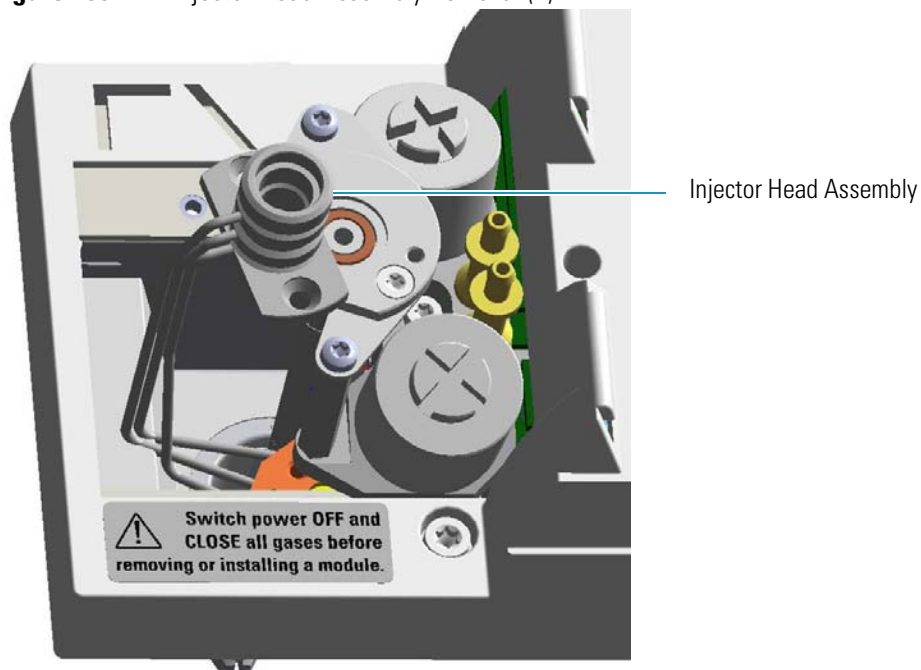

- c. Extract the assembly sideways from its housing. See [Figure 239](#).

**Note** Do not remove the carrier, split, and purge lines' o-rings on the plate.

**Figure 239.** PTVBKF Injector Head Assembly Removal (2)

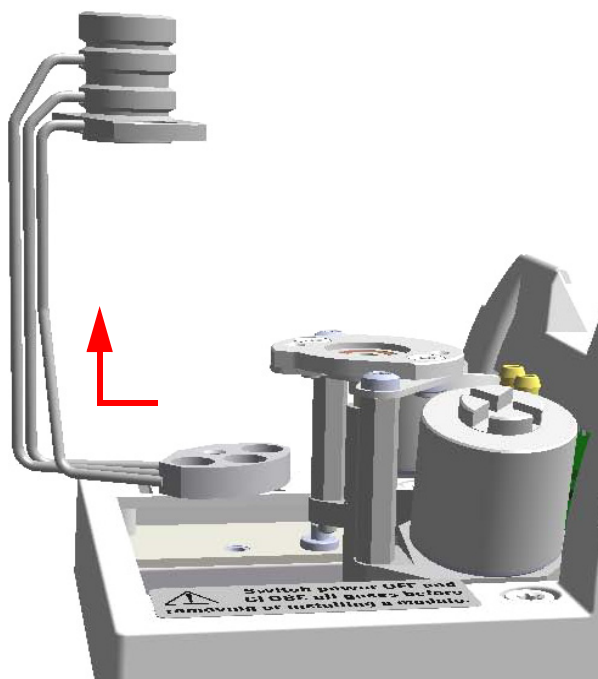

12. Clean the injector head.
  - a. Put the injector head assembly into an ultrasonic cleaner filled with a methanol/acetone mixture (1:1).
  - b. Sonicate the injector head assembly for about half an hour.
  - c. Using tweezers, remove the injector head assembly from the bath and dry it with compressed clean air.
13. Reinstall the injector head.
  - a. Carefully reinstall the injector head assembly into its housing proceeding in the reverse order in which it was removed.
  - b. Screw the two fixing screws and the captive screw fixing the top, and the bottom of the injector head assembly respectively.
14. Remount the fan assembly.
  - a. Carefully reinsert the fan assembly in its seat into the module.
  - b. Move the assembly slightly to position the slot on the right base of the support under the head of the right fixing screw, and for aligning the hole on the left base with the corresponding fixing hole.
  - c. Using tweezers, reinsert the left fixing screw previously removed.
  - d. Using a T10 Torxhead screwdriver, tighten both the left and right fixing screws.
15. Reinsert the liner.
  - a. Place the liner seal over the liner.
  - b. Using tweezers, insert the liner into the injector and push it gently towards the bottom of the injector.
  - c. Reinsert and screw the liner cap using the slotted stubby driver.
16. Reinstall the upper parts of the injector.
  - a. Place the septum into the injector head assembly.



**CAUTION** We suggest replacing the septum with a new one.

- b. Screw and tighten the septum cap to fingertight.



**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

17. Close the module flap cover.
18. If present, move the autosampler towards the module to restore the original alignment.



19. Open the gas supplies.
20. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
21. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
22. Set the normal oven, injectors, and detectors working conditions.

## **6 Injectors Advanced Maintenance**

Cleaning the PTVBKF Injector Head Assembly

## Detectors Advanced Maintenance

This chapter describes TRACE 1300/TRACE 1310 detector modules that do not require routine maintenance, but troubleshooting may indicate they need to be cleaned or replaced.


### Contents

- Removing/Replacing a Detector Module
- Measuring the FID Gas Flows
- Cleaning or Replacing the FID Jet
- Measuring the NPD Gas Flows
- Measuring the FPD Gas Flows
- Cleaning or Replacing the NPD Jet
- Cleaning or Replacing the FPD Mirror Metal Plug
- Cleaning or Replacing the FPD Filter-side Heat Shields
- Cleaning or Replacing the FPD Flame-side Heat Shields
- Replacing the FPD Photomultiplier Tube


## Removing/Replacing a Detector Module

This section provides the instruction for removing/replacing a front/back detector module.


**Figure 240.** Example of Detector Module

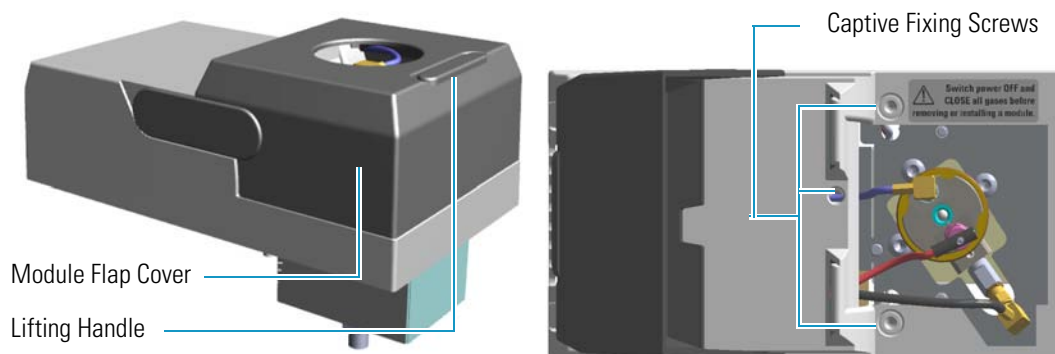


**Figure 241.** Replacing a Front Detector Module

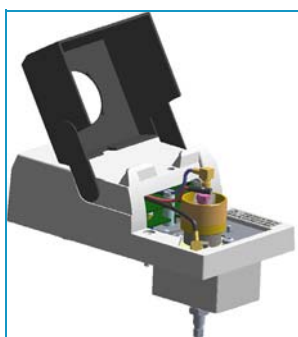


Example of Detector Module

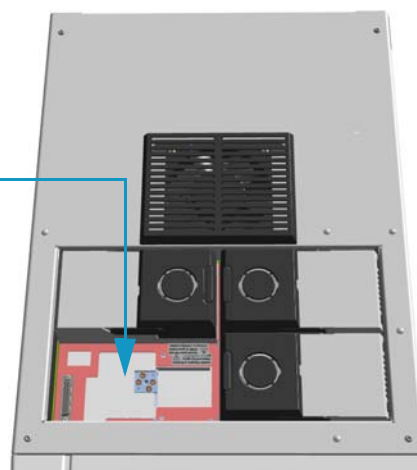

**Figure 242.** Replacing a Back Detector Module

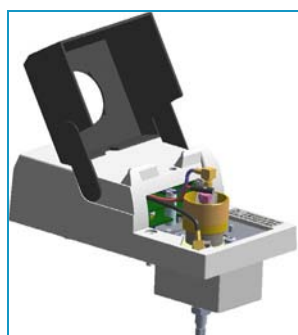

Example of Detector Module

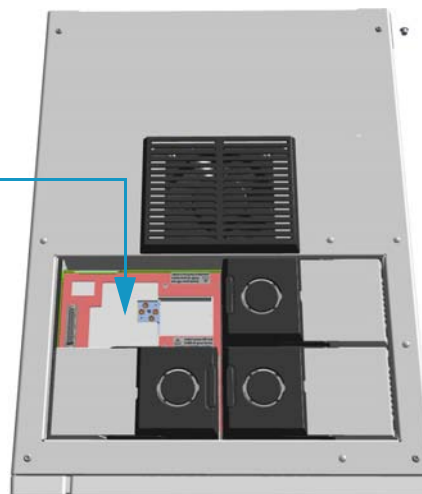


❖ **To remove/replace a front/back detector module**

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.
3. Close the gas supplies.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector on the back of the GC, and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
6. Remove the column from the detector.
  - a. Open the front door of the GC.
  - b. Loosen the retaining nut from the detector fitting on the upper interior wall of the GC oven.
  - c. Remove the analytical column with its nut and ferrule from the bottom of the detector.
7. Remove the detector module.
  - a. Open the module flap cover.
  - b. Using a T20 Torxhead screwdriver, unscrew the three captive fixing screws.
  - c. Keeping the module flap cover open, lift up the module from its seat in the injector housing.
  - d. Place the module on a clean surface.

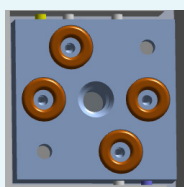

**WARNING** Make sure all the four o-rings are placed into their seats on the gas connection.

**Do not install the module if the o-rings are missing.**

8. Replace the detector module.
  - a. Open the module flap cover.

- b. Keeping the module flap cover open, place it in its seat. Be sure to insert the 25-pin male connector, on the bottom of the module, into the 25-pin female connector on the injector seat of the injector housing.
- c. Use a T20 Torxhead screwdriver to tighten the three captive fixing screws without overtightening.



**CAUTION** To maintain the correct alignment the screws must be tightened in turn. Tighten each screw only a small amount before moving to the next screw. Repeat until all are secure.

- d. Close the module flap cover.
9. Remove the plug, and reconnect the column end to the detector.
10. Open the gas supplies.
11. Check the gas supply for leaks.
  - a. Use a handheld electronic leak detector (Thermo Scientific GLD Pro leak detector or equivalent) to check each fitting for leaks.
  - b. If you detect a leak, tighten the connection and retest it.
  - c. Repeat this process until all connections are leak free.
12. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
13. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
14. Check the module gas connections for leaks.
15. Close the front door of the GC.

## Measuring the FID Gas Flows

Use this procedure for manually measuring the FID fuel and make-up gases.

### ❖ To manually measure the FID flow rates

#### Materials needed

Forceps or tweezers

T10 Torxhead screwdriver

Measuring Tool - Flowmeter Adapter


Electronic flowmeter (Thermo Scientific GFM Pro Flowmeter, or equivalent)

1. Put the GC in standby condition.
2. Turn the flame off. The fuel gases, hydrogen, and air, are automatically closed.
3. Cool the detector to room temperature.
4. Turn the make-up gas off.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

5. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC, and from the wall outlet.
6. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet
7. Open the module flap cover.
8. Disconnect the signal, glow-plug, and polarizing cables from their contacts on the cell top cover. See [Figure 243](#) and [Figure 244](#).

**Figure 243.** FID Cables

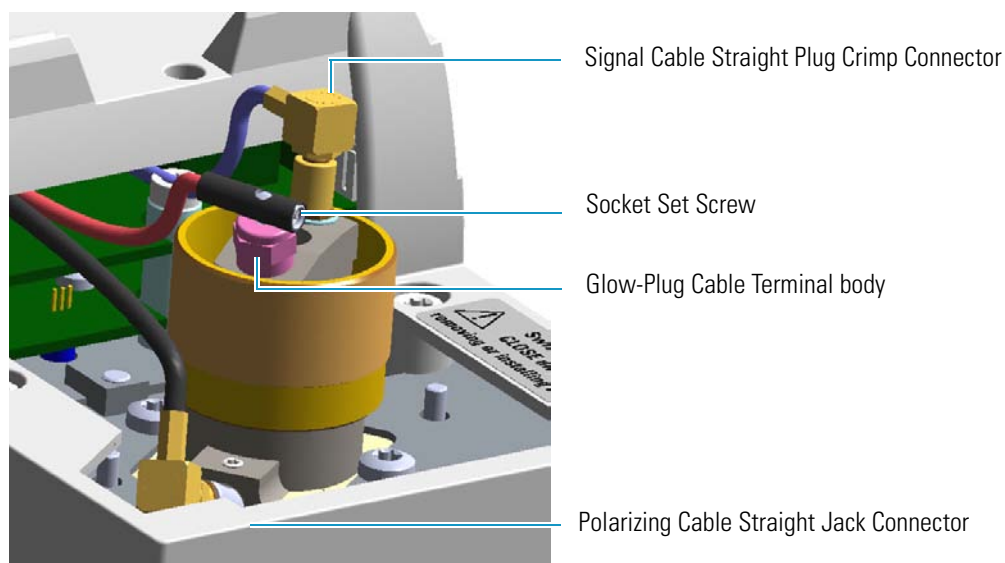


**Figure 244.** FID Cables Disconnection

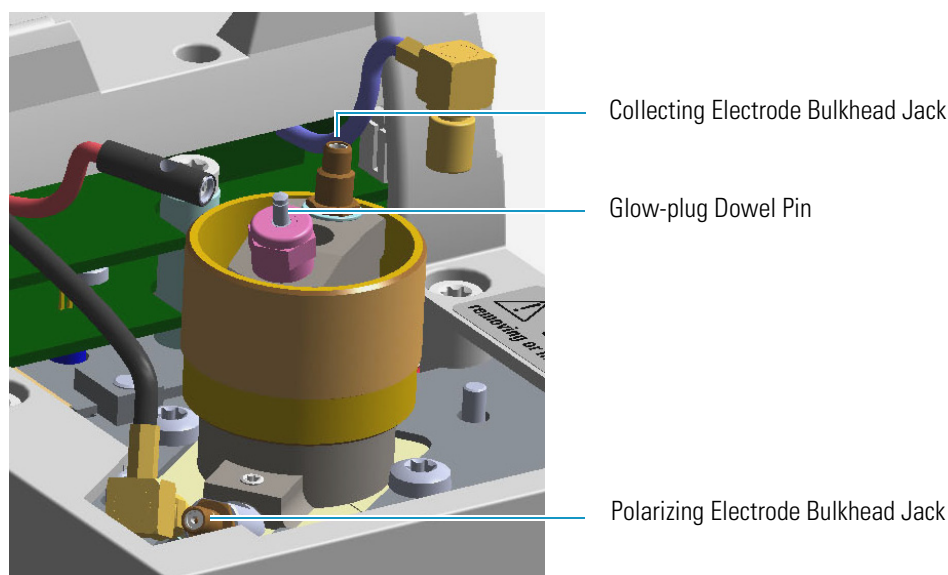

- a. Loosen the glow-plug cable socket set screw using a T10 Torxhead screwdriver. Carefully pull out the terminal body of the glow-plug cable from the dowel pin on the glow-plug element.
  - b. Unscrew and pull out the straight plug crimp connector of the signal cable from the collecting electrode bulkhead jack.
  - c. Pull out the quick coupling straight jack connector of the polarizing cable from the polarizing electrode bulkhead jack.
  - d. Carefully move the cables in order to have free space for handling the detector.
9. Remove the top parts of the detector. See [Figure 245](#).



**Figure 245.** FID Top Parts Removal

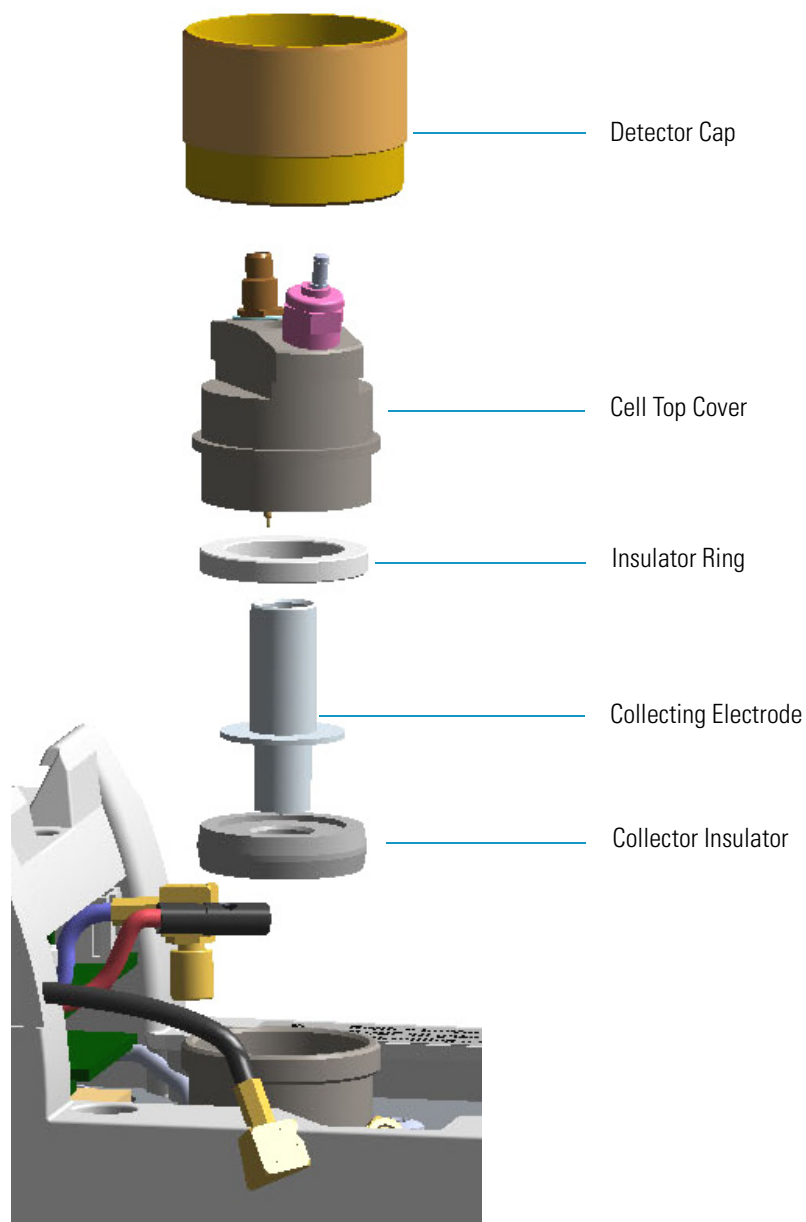


- a. Unscrew and remove the detector cap, paying attention to not rotate the cell top cover.
  - b. Remove the cell top cover and put it in a safe place. Do not damage the FID collecting electrode pin.
10. Remove the electrode. See [Figure 245](#).
- a. Using forceps or tweezers, carefully extract the collecting electrode, the insulator ring and the collector insulator through the top of the detector cell. Place all the parts removed on a clean surface.



**WARNING** The collector insulator must be stored somewhere it can stay clean.


- b. Loosen the polarizing electrode fixing screw using a T10 Torxhead screwdriver, and extract the polarizing electrode from its seat.
  - c. Place the polarizing electrode on a clean surface.
11. Insert the FID gases measuring tool. See [Figure 246](#).
- a. Insert the measuring tool into the detector body.
  - b. Push the measuring tool downwards until it stops.

**Figure 246.** FID Gases Measuring Tool

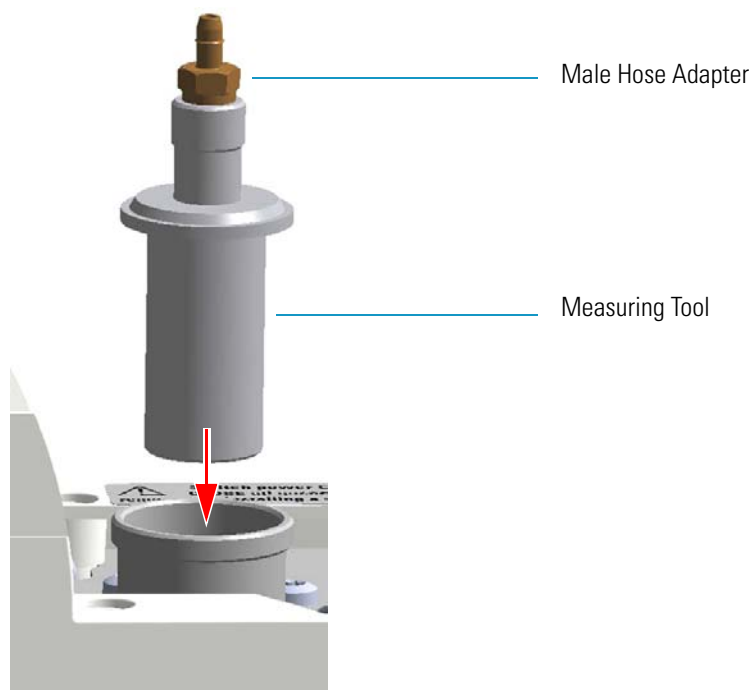


- c. Screw the detector cap onto the measuring tool. See [Figure 247](#).

**Figure 247.** FID Gases Measuring Tool Assembling

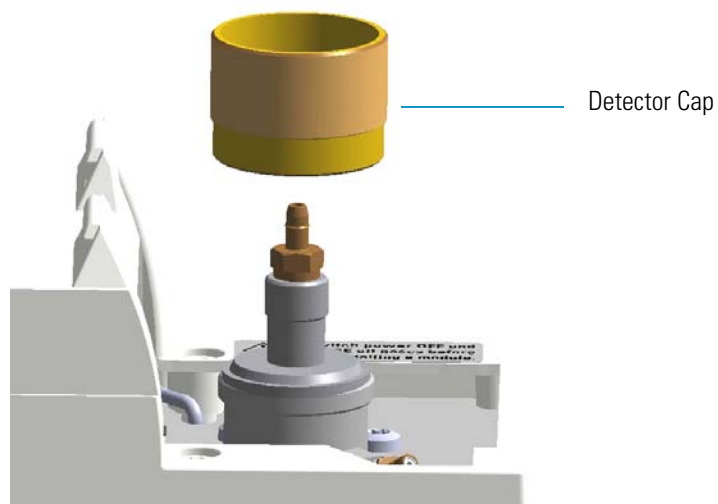

- d. Place the flowmeter tubing over the male hose adapter of the measuring tool to begin measuring flows. See [Figure 248](#).

**Figure 248.** Measuring Tool to Flowmeter Connection

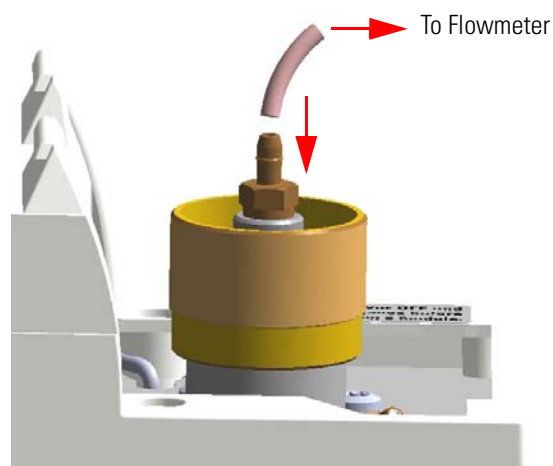

12. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
13. Measure the gas flow rate.
  - a. Open the GC user interface.
  - b. Turn the **Hydrogen** flow On.
    - i. Measure the gas flow and verify that the actual flow rate corresponds to the setpoint value.

- c. Turn the hydrogen flow Off.
  - d. Turn the **Air** supply On.
    - i. Measure the gas flow and verify that the actual flow rate corresponds to the setpoint value.
  - e. Turn the Air flow Off.
  - f. Turn the **Makeup** gas flow On.
    - i. Measure the gas flow and verify that the actual flow rate corresponds to the setpoint value.
  - g. Turn the Makeup Off.
14. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC, and from the wall outlet.
15. Unscrew the detector cap, then remove the measuring tool from the detector body.
16. Reinstall the electrodes.
  - a. Reinsert the polarizing electrode into its seat pushing while screwing it, then screw the fixing screw.
  - b. Remount the collector insulator, the collecting electrode, and its insulator ring.
17. Remount the top parts of the detector.
  - a. Remount the cell top cover.
  - b. Screw the detector cap paying attention to not rotate the cell top cover.
  - c. Reconnect the signal, glow-plug, and polarizing cables.
18. Close the module flap cover.
19. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
20. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
21. Turn the make-up gas on.
22. Set the normal oven and detector working conditions.
23. Ignite the flame. The fuel gases, hydrogen and air, are automatically opened.

## Cleaning or Replacing the FID Jet

It is good practice to clean the jet semiannually when you are analyzing particularly dirty compounds such that the flame does not burn properly. In case of actual damage, the jet must be replaced.

### ❖ To clean or replace the FID ceramic jet

#### Materials needed

GC-grade methanol

Distilled water

Screwdriver?

Forceps or tweezers

T10 Torxhead screwdriver

Elbowed box wrench, 8-mm


Ceramic jet (if necessary)

1. Put the GC in standby condition.
2. Turn the flame off. The fuel gases, hydrogen, and air, are automatically closed.
3. Cool the detector to room temperature.
4. Turn the make-up gas off.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

5. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC, and from the wall outlet.
  - c. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
6. Open the module flap cover.
7. Disconnect the signal, glow-plug, and polarizing cables from their contacts on the cell top cover. See [Figure 249](#) and [Figure 250](#).

**Figure 249.** FID Cables

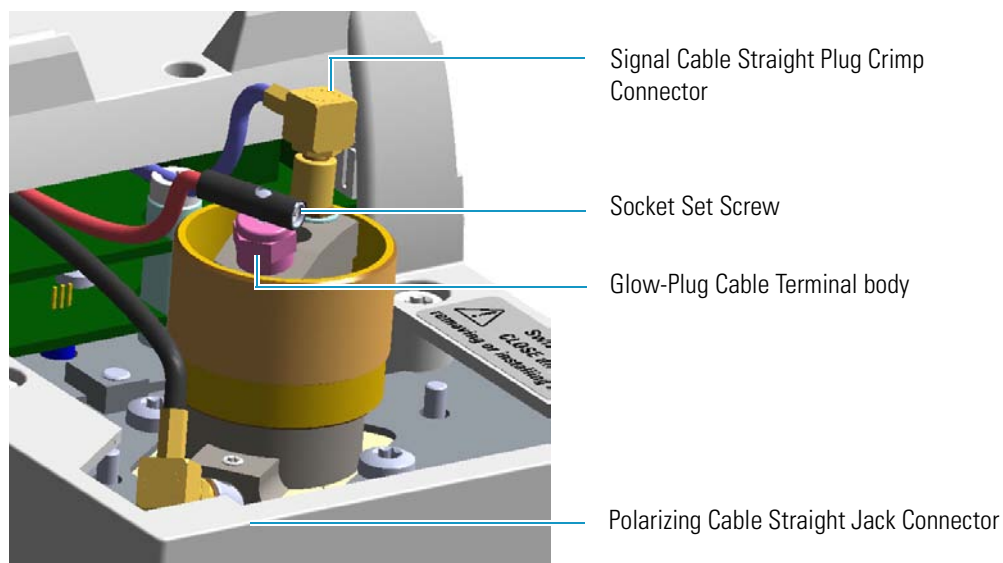


**Figure 250.** FID Cables Disconnection

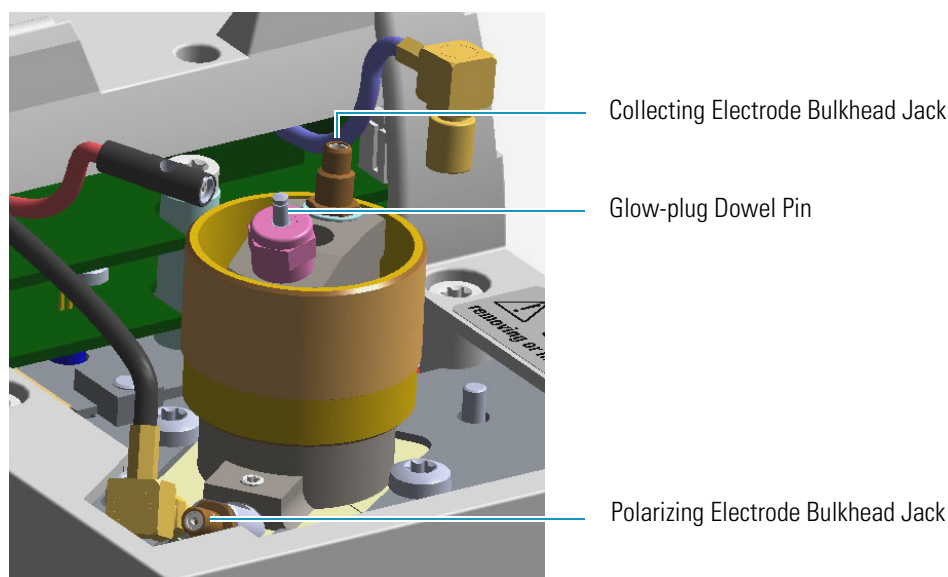

- a. Loosen the glow-plug cable socket set screw using a T10 Torxhead screwdriver. Carefully pull out the terminal body of the glow-plug cable from the dowel pin on the glow-plug element.
  - b. Unscrew and pull out the coaxial straight plug crimp connector of the signal cable from the collecting electrode bulkhead jack.
  - c. Pull out the quick coupling coaxial straight jack connector of the polarizing cable from the polarizing electrode bulkhead jack.
  - d. Carefully move the cables in order to have free space for handling the detector.
8. Remove the top parts of the detector. See [Figure 251](#).

**Figure 251.** FID Top Parts Removal

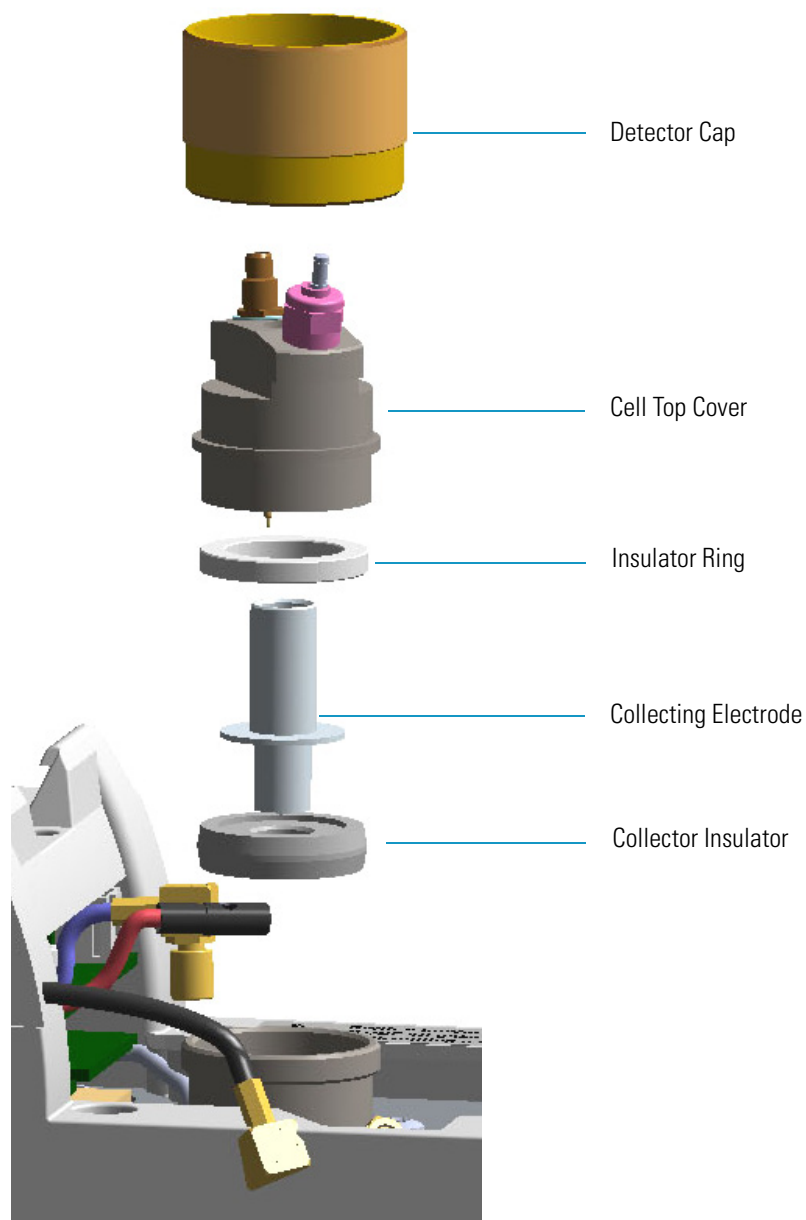


- a. Unscrew and remove the detector cap paying attention to not rotate the cell top cover.
  - b. Remove the cell top cover and put it in a safe place. Do not damage the FID collecting electrode pin.
9. Remove the electrode. See [Figure 251](#).
- a. Using forceps or tweezers, carefully extract the collecting electrode, the insulator ring, and the collector insulator through the top of the detector cell. Place all the parts removed on a clean surface.



**WARNING** The collector insulator must be stored somewhere it can stay clean.


- b. Loosen the polarizing electrode fixing screw using a T10 Torxhead screwdriver, and extract the polarizing electrode from its seat. See [Figure 252](#).

**Figure 252.** FID Jet (2)

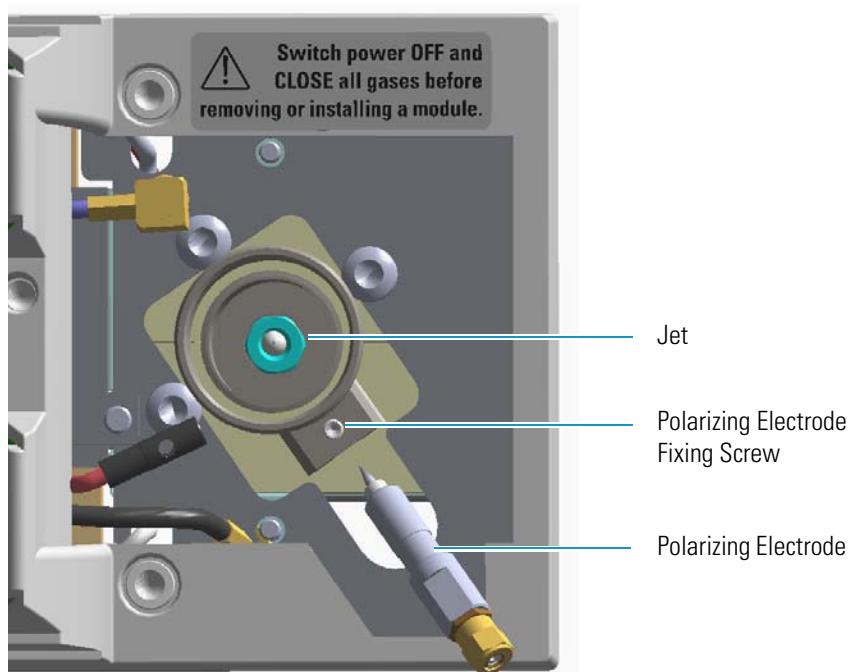

- c. Place the polarizing electrode on a clean surface.
10. Remove, clean, and reinstall the jet.
    - a. Use the elbowed box wrench for unscrewing and extracting the jet through the detector body. See [Figure 253](#).

**Figure 253.** FID Jet (2)

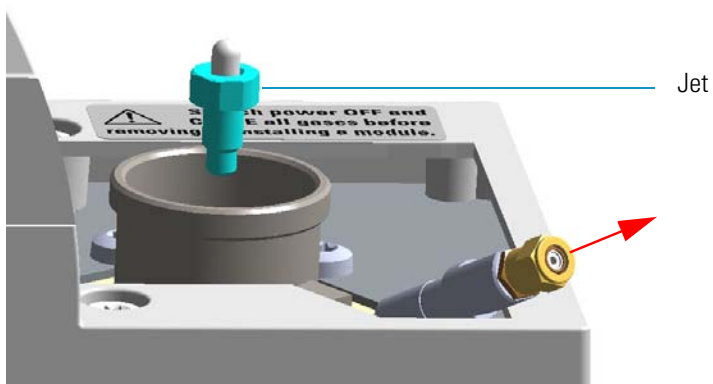



- b. Place the jet in the ultrasonic cleaner, filled with liquid detergent.
- c. Sonicate the jet for about five minutes.
- d. Handling the jet with forceps or tweezers. Rinse the jet using distilled water, then methanol.
- e. Place the jet on a paper towel, and let the jet air dry.

**Note** If after the cleaning physical condition of the jet does not permit its reuse, replace the jet with a new one.

- f. Using the proper tool, reinsert and screw the cleaned (or a new) jet into its seat in the detector body.
11. Reinstall the electrodes.
- a. Reinsert the polarizing electrode into its seat pushing while screwing it, then screw the fixing screw.
  - b. Remount the collector insulator, the collecting electrode, and its insulator ring.
12. Remount the top parts of the detector.
- a. Remount the cell top cover.
  - b. Screw the detector cap paying attention to not rotate the cell top cover.
  - c. Reconnect signal, glow-plug, and polarizing cables.
13. Close the module flap cover.
14. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
15. Power on the GC.
- a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
16. Turn the make-up gas on.
17. Set the normal detector working conditions.
18. Ignite the flame. The fuel gases, hydrogen, and air, are automatically opened.

## Measuring the NPD Gas Flows

Use this procedure for manually measuring the NPD fuel and make-up gases.

### ❖ To manually measure the NPD flow rates

#### Materials needed

Forceps or tweezers

T6 Torxhead key

T10 Torxhead screwdriver

6 mm wrench

Measuring Tool Flowmeter Adapter

Electronic flowmeter (Thermo Scientific GFM Pro Flowmeter, or equivalent)

1. Put the GC in standby condition.
2. Switch off the thermionic source. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.
4. Turn the make-up gas off.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

5. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
6. Unplug the power cable from the AC input connector of the **NPD Thermionic Source Power Module** and from the wall outlet.

If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.


7. Unplug the signal cable from the detector module. See [Figure 254](#).

**Figure 254.** NPD Signal Cable Removal

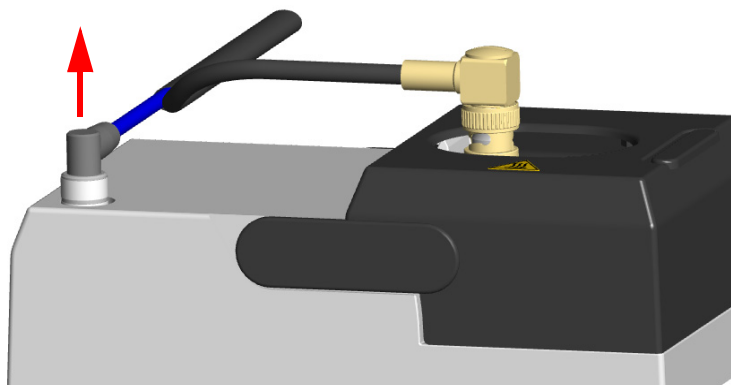


8. Open the module flap cover.
9. Remove the thermionic source assembly cable. See [Figure 255](#).

**Figure 255.** NPD Thermionic Source Cable Removal

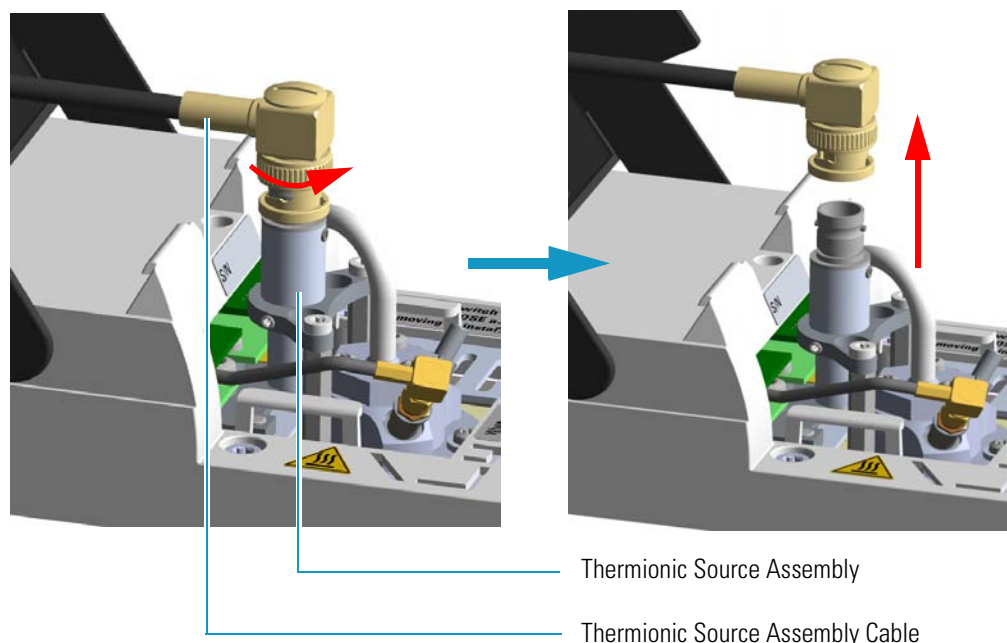

- a. Twist the ring to disconnect the thermionic source assembly cable.
  - b. Push and twist the lock so that the button slides up in the groove, then pull the cable ends apart.
10. Unscrew and pull out the straight jack connector of the collecting electrode cable from the collecting electrode bulkhead jack. See [Figure 256](#).

**Figure 256.** NPD Collecting Electrode Cable Removal

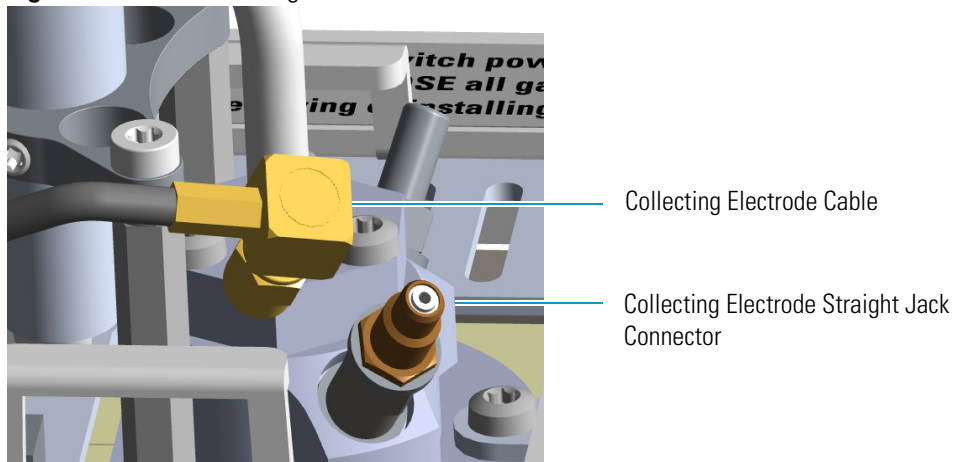


11. Remove the thermionic source



**CAUTION** The thermionic source is delicate. Be careful not to break or crack the source. When performing maintenance on the NPD, avoid touching the source with your fingers, and prevent it from coming in contact with other surfaces. Wear clean, lint-free gloves to prevent contamination of parts with dirt and skin oils.


- a. Using a T6 Torxhead key, loosen the dowel which fix the thermionic source connector to the thermionic source support. See [Figure 257](#).

**Figure 257.** NPD Thermionic Source Removal (1)

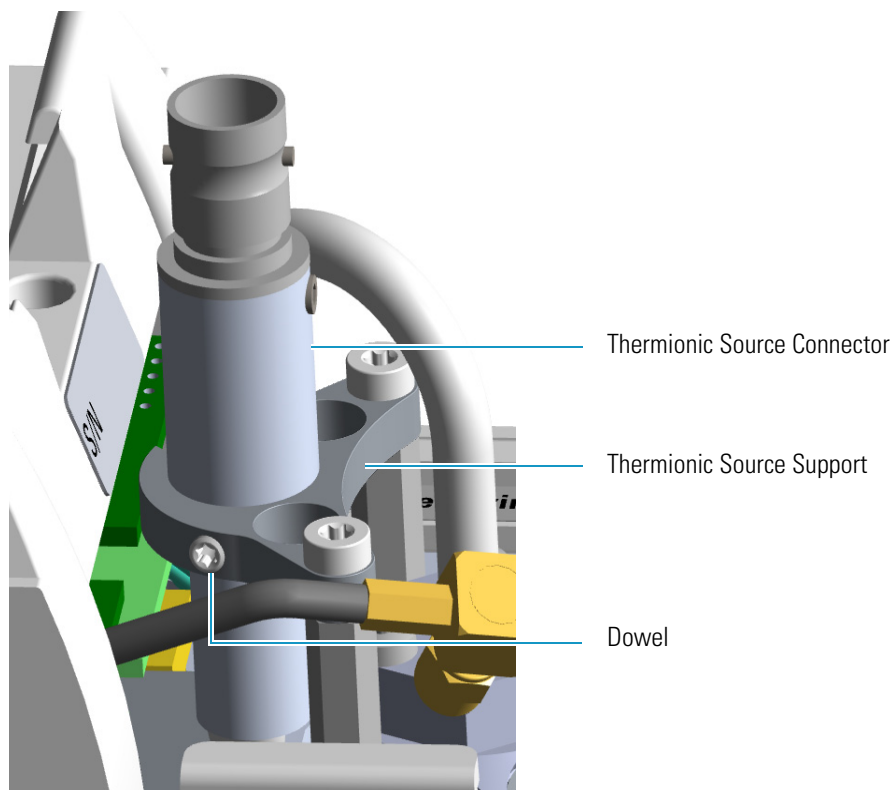


- b. Using the T10 Torxhead screwdriver, remove the three T10 Torx screws from the thermionic source assembly. See [Figure 258](#).

**Figure 258.** NPD Thermionic Source Removal (2)

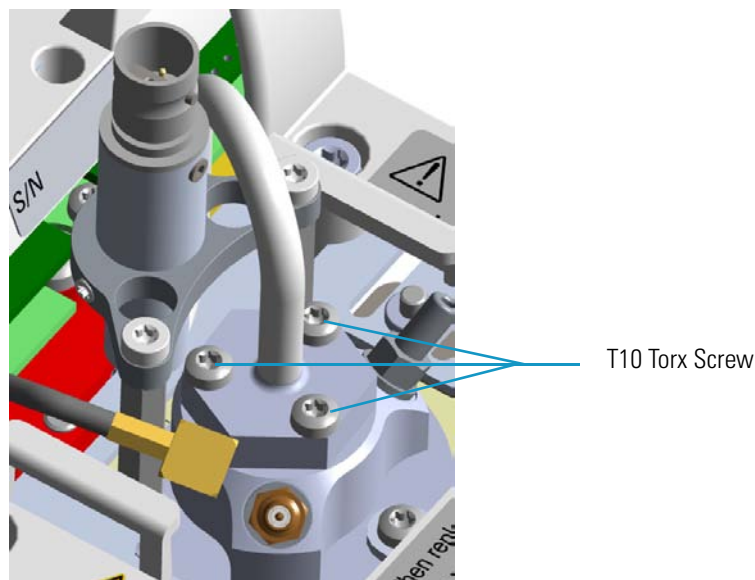

- c. Gently lift up the thermionic source connector guiding the flexible cable from the thermionic source support, then remove the thermionic source assembly from the NPD body. Avoid bumping the bead on the sides of the collector. See [Figure 259](#).


**Figure 259.** NPD Thermionic Source Removal (3)

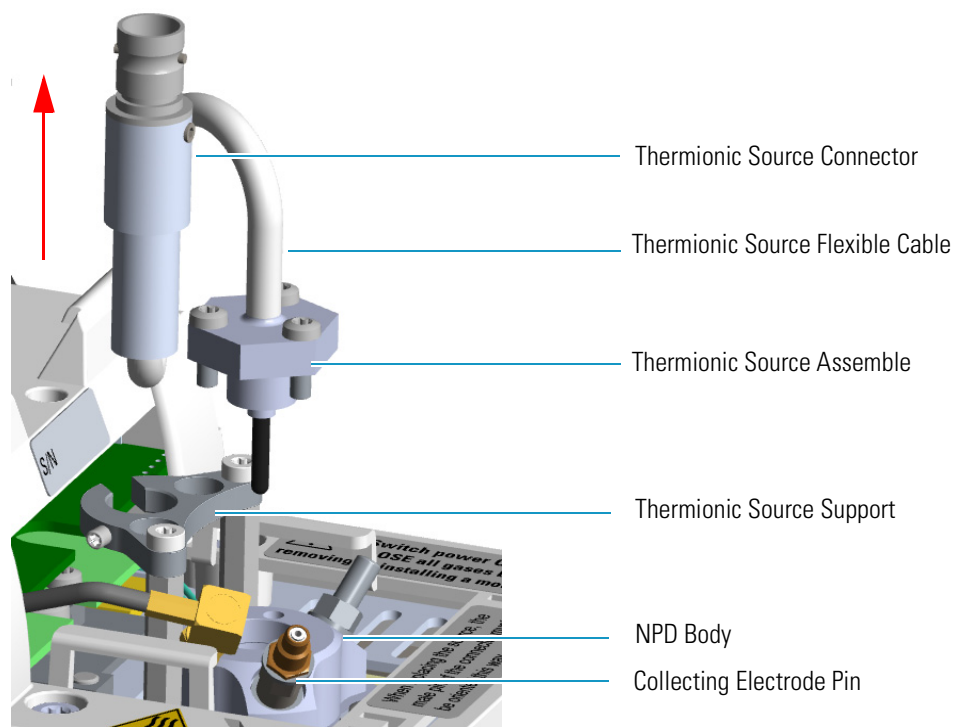

12. Remove the collecting electrode and the insulator. See [Figure 260](#).

- Using a 6-mm wrench, unscrew and remove the collecting electrode pin from its seat on the detector body.
- Using forceps or tweezers, extract the collecting electrode through the top of the detector body.
- Using forceps or tweezers, extract the insulator through the top of the detector body. Place the insulator on a clean surface.


**Figure 260.** NPD Insulator and Collecting Electrode Removal

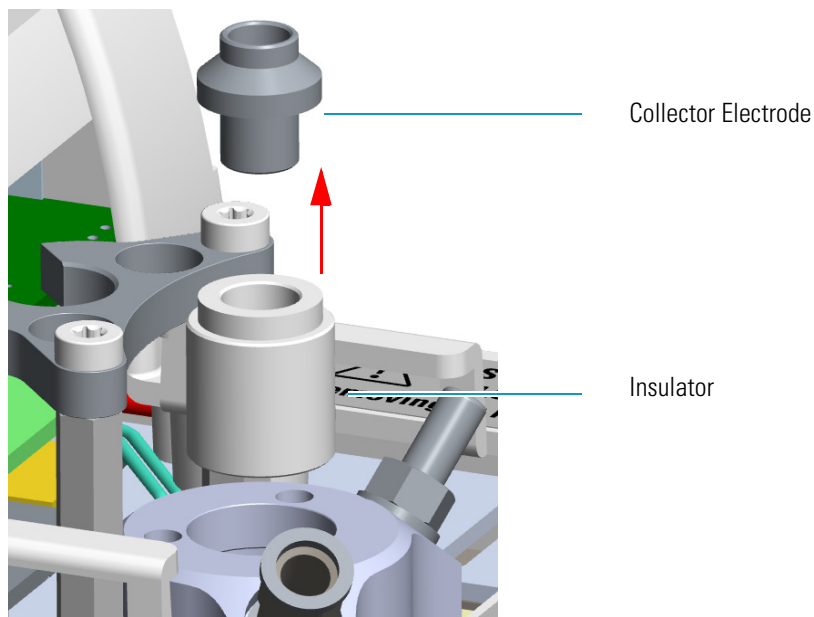

**WARNING** The collector insulator must be stored somewhere it can stay clean.


13. Insert the NPD gases measuring tool. See [Figure 261](#).

**Figure 261.** NPD Gases Measuring Tool (1)

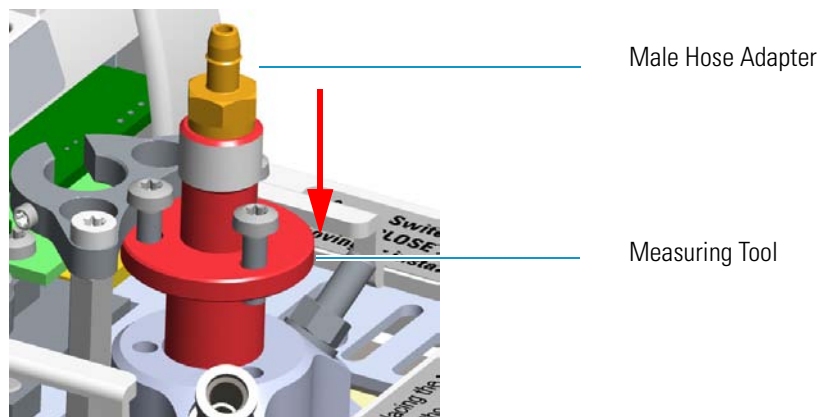

- a. Insert the flowmeter adapter measuring tool into the detector body.
- b. Push the flowmeter adapter measuring tool downwards until it stops.
- c. Rotate the measuring tool to align the fixing holes, then insert the three screws used to fix the thermionic source.
- d. Fix the flowmeter adapter measuring tool to the detector body using a T10 Torxhead screwdriver.
- e. Place the flowmeter tubing over the male hose adapter of the measuring tool to begin measuring flows. See [Figure 262](#).


**Figure 262.** NPD Gases Measuring Tool (2)

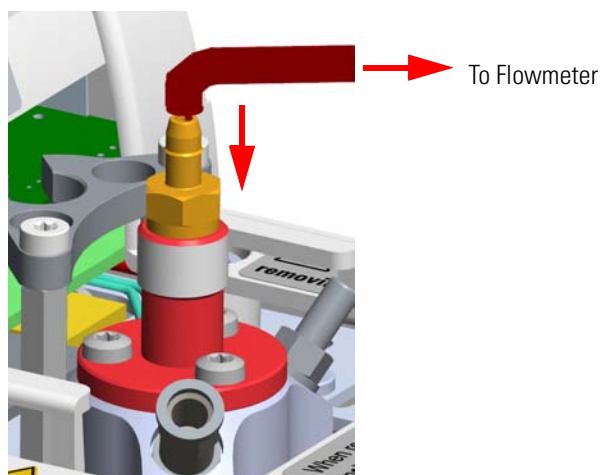

14. Measure the gas flow rate.
  - a. Open the GC user interface.
  - b. Turn the **Hydrogen** flow On.
    - i. Measure the gas flow and verify that the actual flow rate corresponds to the setpoint value.
  - c. Turn the hydrogen flow Off.
  - d. Turn the **Air** supply On.
    - i. Measure the gas flow and verify that the actual flow rate corresponds to the setpoint value.
  - e. Turn the Air flow Off.
  - f. Turn the **Makeup** gas flow On.
    - i. Measure the gas flow and verify that the actual flow rate corresponds to the setpoint value.
  - g. Turn the Makeup Off.
15. Remove the flowmeter adapter measuring tool from the detector body proceeding in the reverse order in which it was installed.

- a. Using forceps or tweezers, reinstall the insulator into the detector body.
16. Reinstall the insulator and the collecting electrode.
  - a. Using forceps or tweezers, reinstall the insulator into the detector body.
  - b. Place the collecting electrode in its housing.
  - c. Using a 6-mm wrench, screw the collecting electrode pin into its seat on the detector body.
  - d. Screw the straight plug crimp connector of the collecting electrode cable to the collecting electrode bulkhead jack.
17. Reinstall the thermionic source.
  - a. Remove the protective cap covering the new thermionic source.
  - b. Mount and guide the new source assembly on the NPD body proceeding in the reverse order in which it was removed. Be careful not to bump the bead on the sides of the body and collecting electrode.
  - c. Rotate and align the thermionic source connector, then tighten the dowel to fix the connector to the support using the T6 Torxhead key. See [Figure 263](#).


**Figure 263.** Reinstall Thermionic Source (1)

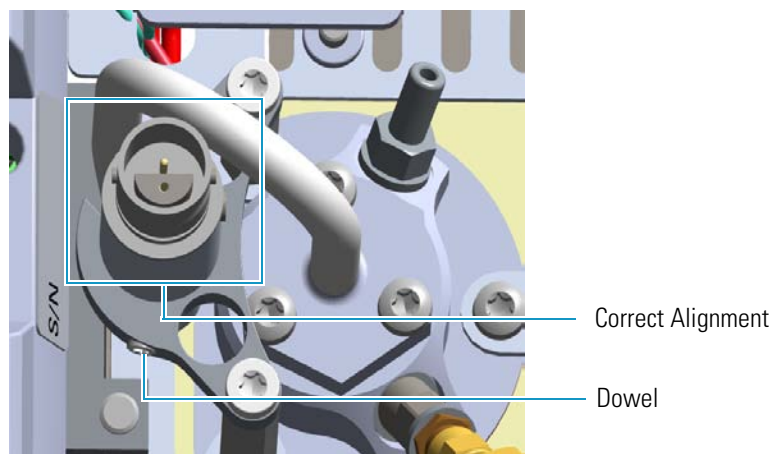

18. Reconnect and screw the straight jack connector of the collecting electrode cable to the collecting electrode bulkhead jack. See [Figure 264](#).




**Figure 264.** Reinstall Thermionic Source (2)

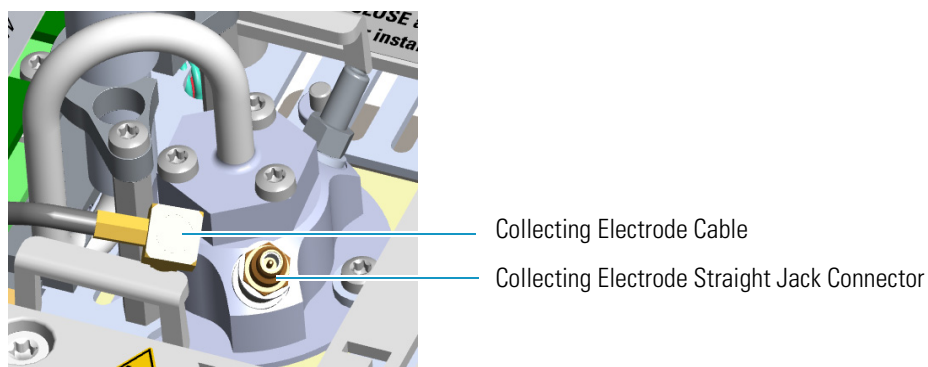


19. Reconnect the source assembly cable to the NPD cable and twist the ring to lock the connection. See [Figure 265](#).

**Figure 265.** Reinstall Thermionic Source (3)

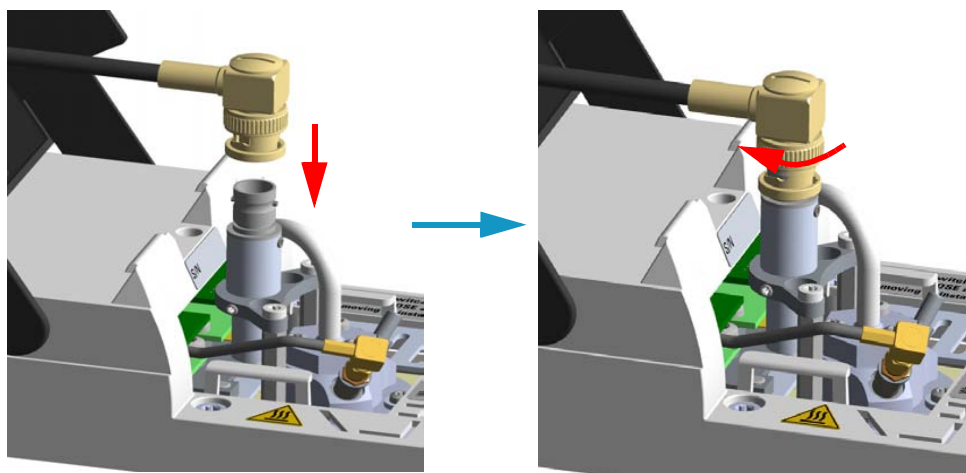

**Figure 266.** Plug NPD Signal Cable

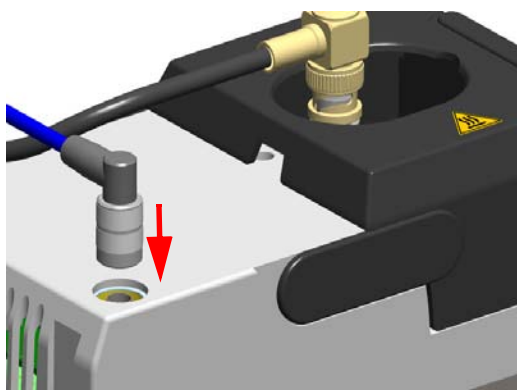

22. Plug the power cable to the AC input connector of the **NPD Thermionic Source Power Module** and to the wall outlet. The LED marked **On** blinks green.

If other external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.

23. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
24. Turn the make-up gas on.
25. Heat the detector the working temperature.
26. Switch on the thermionic source. Hydrogen and air are automatically opened. The LED marked **On** of the NPD Thermionic Source Power Module becomes solid green.

## Cleaning or Replacing the NPD Jet

You should cleaning the jet semiannually when you analyze particularly dirty compounds that the flame does not burnt properly. In case of breakage or cracking, the jet must be replaced.

### ❖ To clean or replace the NPD jet

#### Materials needed

|                          |
|--------------------------|
| Ultrasonic cleaner       |
| Liquid detergent         |
| GC-grade methanol        |
| Distilled water          |
| Forceps or tweezers      |
| T6 Torxhead key          |
| T10 Torxhead screwdriver |
| 6 mm wrench              |
| Elbowed box wrench, 8-mm |
| Jet (if necessary)       |

1. Put the GC in standby condition.
2. Switch off the thermionic source. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.
4. Turn the make-up gas off.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

5. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
6. Unplug the power cable from the AC input connector of the **NPD Thermionic Source Power Module** and from the wall outlet.

If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.


7. If present, move the autosampler away from the module.
8. Unplug the signal cable from the detector module. See [Figure 267](#).

**Figure 267.** NPD Signal Cable Removal

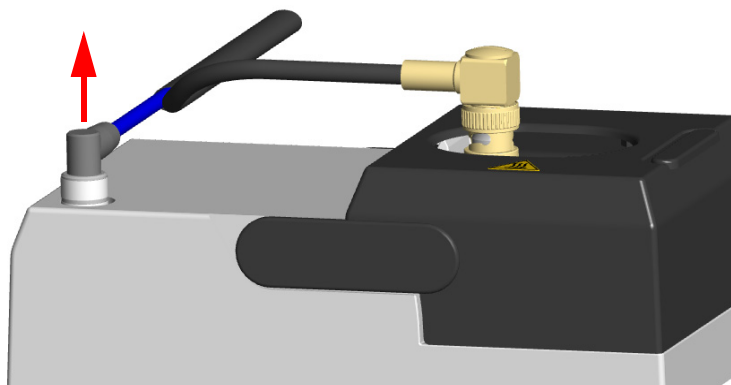


9. Open the module flap cover.
10. Remove the thermionic source assembly cable. See [Figure 268](#).

**Figure 268.** NPD Thermionic Source Cable Removal

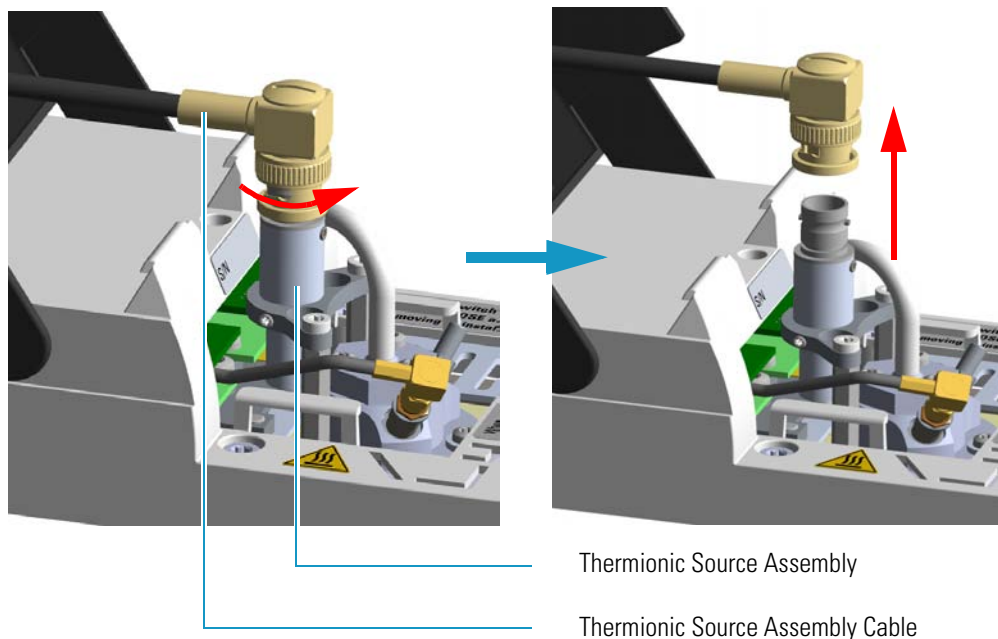

- a. Twist the ring to disconnect the thermionic source assembly cable.
  - b. Push and twist the lock so that the button slides up in the groove, then pull the cable ends apart.
11. Unscrew and pull out the straight jack connector of the collecting electrode cable from the collecting electrode bulkhead jack. See [Figure 269](#).

**Figure 269.** NPD Collecting Electrode Cable Removal

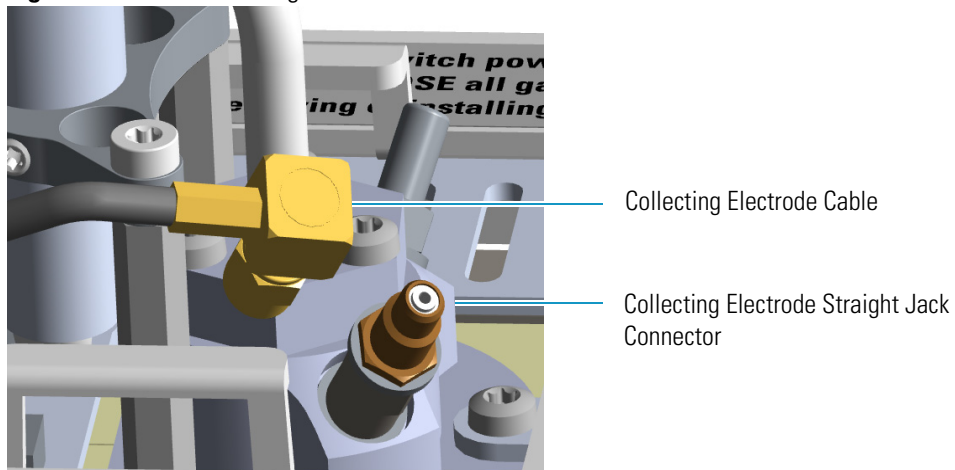


12. Remove the thermionic source



**CAUTION** The thermionic source is delicate. Be careful not to break or crack the source. When performing maintenance on the NPD, avoid touching the source with your fingers, and prevent it from coming in contact with other surfaces. Wear clean, lint-free gloves to prevent contamination of parts with dirt and skin oils.


- a. Using a T6 Torxhead key, loosen the dowel which fixes the thermionic source connector to the thermionic source support. See [Figure 270](#).

**Figure 270.** NPD Thermionic Source Removal (1)

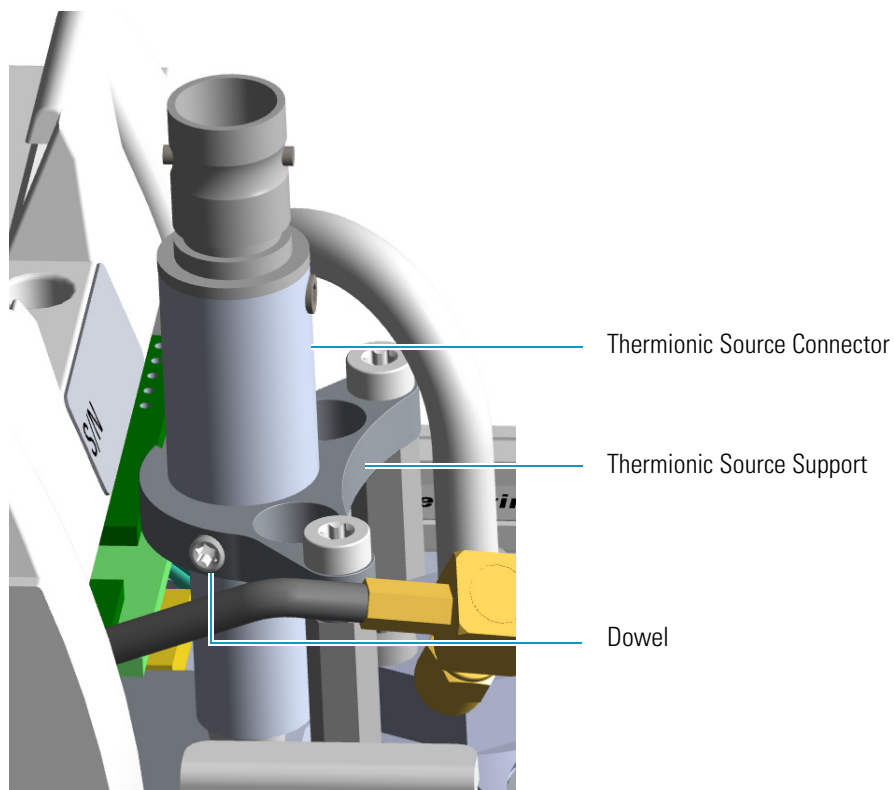


- b. Using the T10 Torxhead screwdriver, remove the three T10 Torx screws from the thermionic source assembly. See [Figure 271](#).

**Figure 271.** NPD Thermionic Source Removal (2)

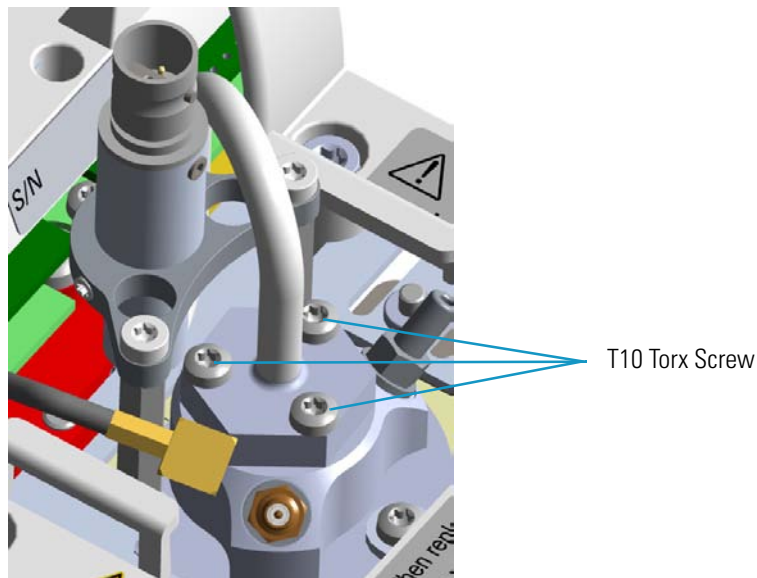


- c. Gently lift up the thermionic source connector guiding the flexible cable from the thermionic source support, then remove the thermionic source assembly from the NPD body. Avoid bumping the bead on the sides of the collector. See [Figure 272](#).

**Figure 272.** NPD Thermionic Source Removal (3)

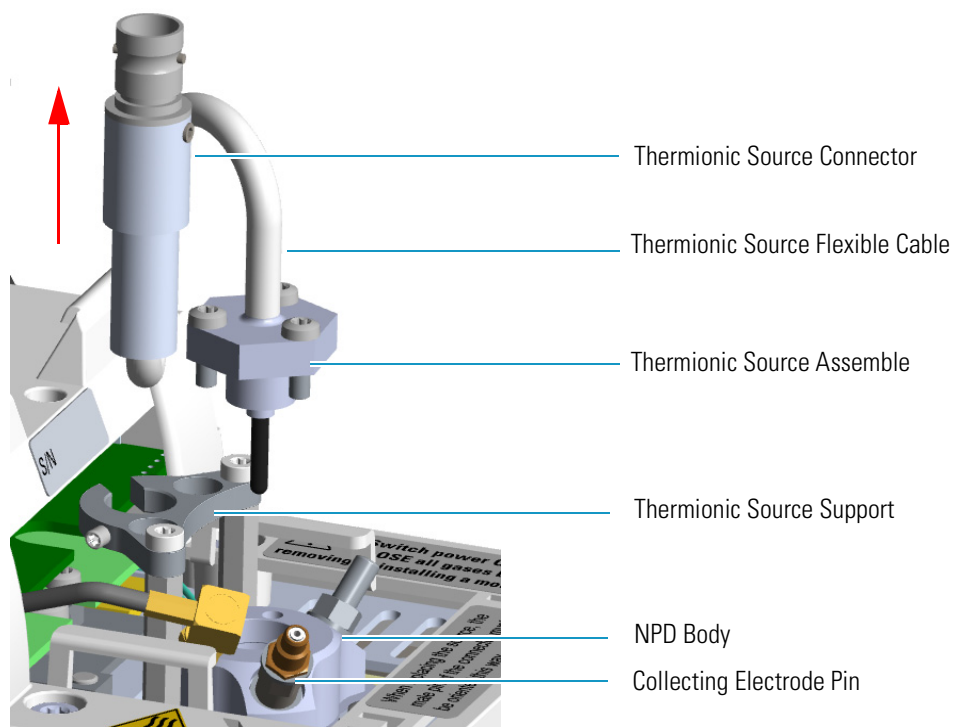

13. Remove the collecting electrode. See [Figure 273](#).
  - a. Using a 6-mm wrench, unscrew and remove the collecting electrode pin from its seat on the detector body.
  - b. Using forceps or tweezers, extract the collecting electrode through the top of the detector body.

**Figure 273.** NPD Collecting Electrode Removal

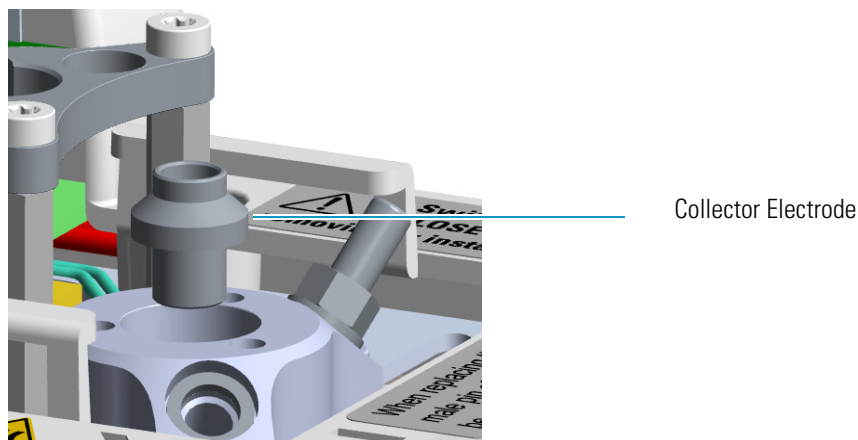


14. Remove and clean the jet. See [Figure 274](#).
  - a. Using forceps or tweezers, extract the insulator through the top of the detector body. Place the insulator on a clean surface.



**WARNING** The collector insulator must be stored somewhere it can stay clean.

- b. Use the elbowed box wrench for unscrewing and extracting the jet through the detector body.

**Figure 274.** Jet Removal

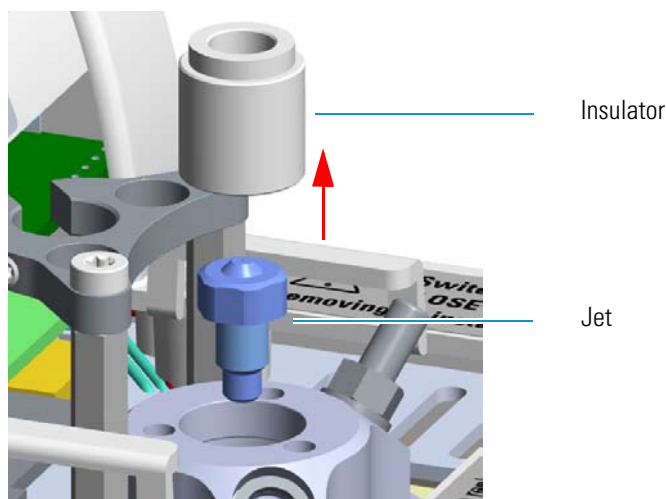

- c. Place the jet in the ultrasonic cleaner, filled with liquid detergent.
- d. Sonicate the jet for about five minutes.

**Note** If after cleaning the physical condition of the jet does not permit its reuse, replace the jet with a new one.

15. Reinsert the jet

- a. Using the proper tool, reinsert the cleaned (or a new) jet, into the detector body and screw it in its seat.
- b. Using forceps or tweezers, reinstall the insulator into the detector body.


16. Reinstall the collecting electrode.

- a. Place the cleaned (or a new) collecting electrode in its housing.
- b. Using a 6-mm wrench, screw the collecting electrode pin into its seat on the detector body.
- c. Screw the straight plug crimp connector of the collecting electrode cable to the collecting electrode bulkhead jack.

17. Reinstall the thermionic source.

- a. Remove the protective cap covering the new thermionic source.
- b. Mount and guide the new source assembly on the NPD body proceeding in the reverse order in which it was removed. Be careful not to bump the bead on the sides of the body and collecting electrode.
- c. Rotate and align the thermionic source connector, then tighten the dowel to fix the connector to the support using the T6 Torxhead key. See [Figure 275](#).


**Figure 275.** Reinstall Thermionic Source (1)

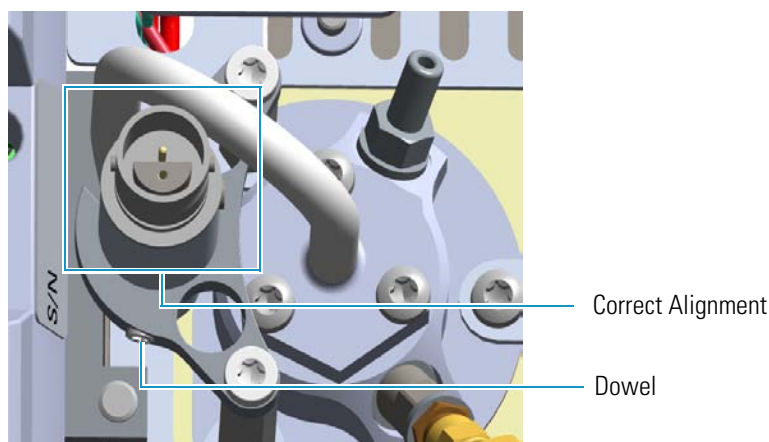

18. Reconnect and screw the straight jack connector of the collecting electrode cable to the collecting electrode bulkhead jack. See [Figure 276](#).




**Figure 276.** Reinstall Thermionic Source (2)

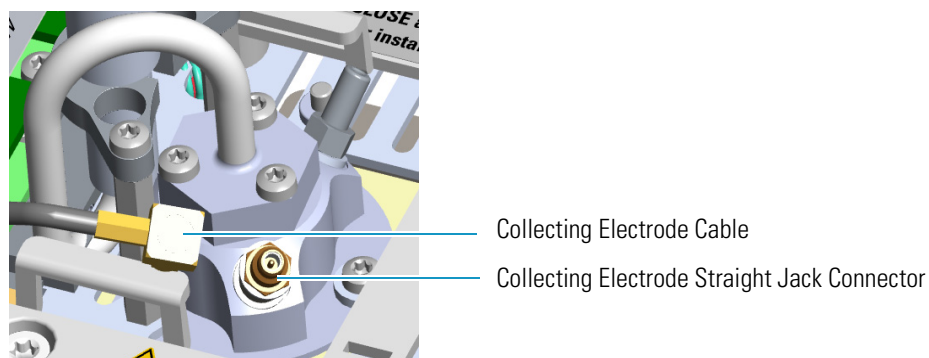


19. Reconnect the source assembly cable to the NPD cable and twist the ring to lock the connection. See [Figure 277](#).

**Figure 277.** Reinstall Thermionic Source (3)

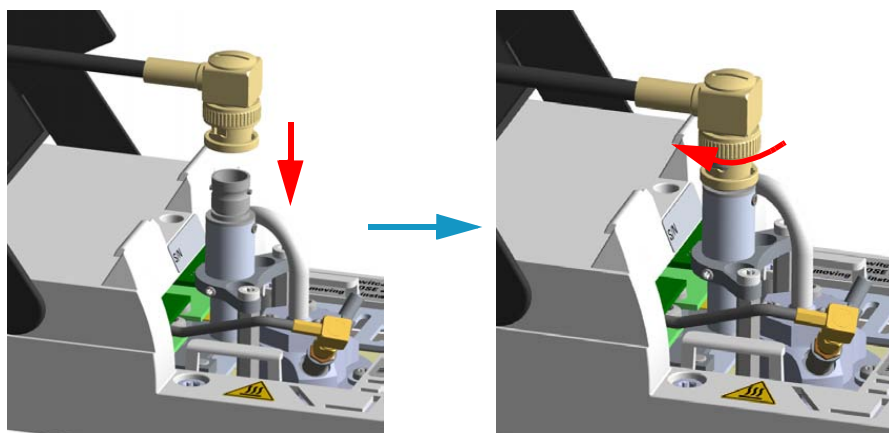

20. Close the module flap cover.
21. Plug in the signal cable into its contact on the detector module. See [Figure 278](#).

**Figure 278.** Plug NPD Signal Cable

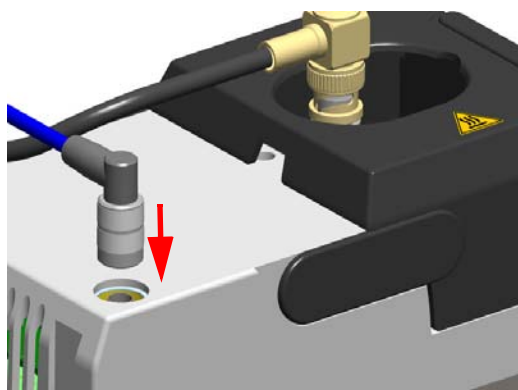

22. If present, move the autosampler towards the module to restore the original alignment.

23. Plug the power cable to the AC input connector of the **NPD Thermionic Source Power Module** and to the wall outlet. The LED marked **On** blinks green.

If other external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.

24. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
25. Turn the make-up gas on.
26. Heat the detector the working temperature.
27. Switch on the thermionic source. Hydrogen and air are automatically opened. The LED marked **On** becomes solid green.

## Measuring the FPD Gas Flows

Use this procedure for manually measuring the FPD fuel gases.

### ❖ To manually measure the FPD flow rates

#### Materials needed

Forceps or tweezers

T10 Torxhead screwdriver

Measuring Tool Flowmeter Adapter


Electronic flowmeter (Thermo Scientific GFM Pro Flowmeter, or equivalent)

1. Put the GC in standby condition.
2. Turn the flame off. The fuel gases, hydrogen, and air, are automatically closed.
3. Cool the detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Remove the column end from the bottom of the detector, then plug the bottom of the detector using the blind cap.
5. Insert the FPD gases measuring tool. See [Figure 279](#).
  - a. Remove the chimney cap.
  - b. Place the measuring tool on the chimney.

**Figure 279.** FPD Gases Measuring Tool

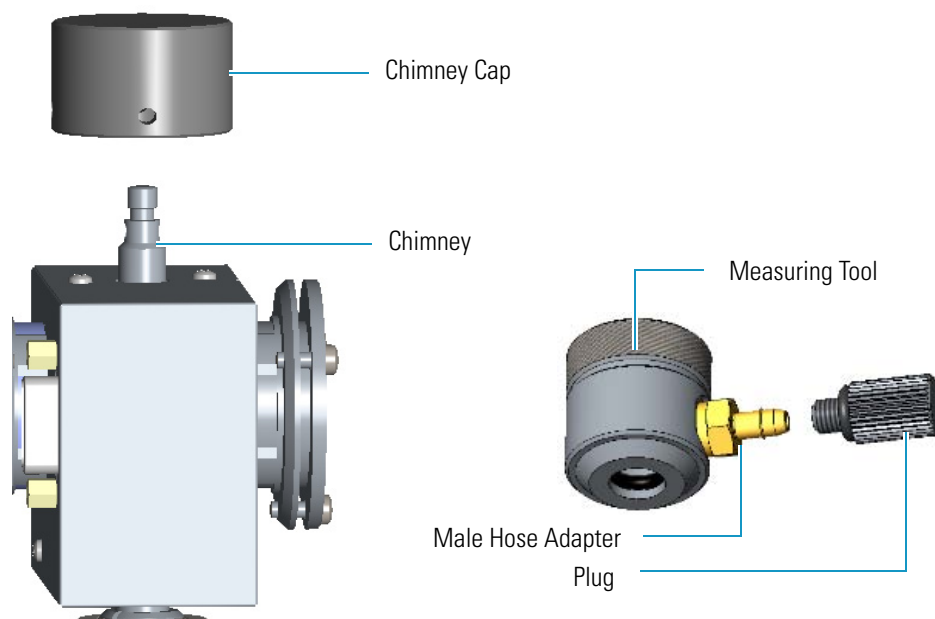

**Note** For measuring the flows use the male hose adapter. The plug must be used instead the male hole adapter when a leak test is required.

6. Place the flowmeter tubing over the male hose adapter of the measuring tool to begin measuring flows.
7. Measure the gas flow rate.
  - a. Open the GC user interface.
  - b. Turn the **Hydrogen** flow On.
    - i. Measure the gas flow and verify that the actual flow rate corresponds to the setpoint value.
  - c. Turn the hydrogen flow Off.
  - d. Turn the **Air** supply On.
    - i. Measure the gas flow and verify that the actual flow rate corresponds to the setpoint value.
  - e. Turn the Air flow Off.
8. Remove the flowmeter tubing from the male hose adapter of the measuring tool.
9. Remove the measuring tool from the chimney.
10. Replace the chimney cap over the chimney.
11. Remove the blind cap from the bottom of the detector.
12. Reconnect the column end into the bottom of the detector.
13. Ignite the flame. The fuel gases, hydrogen and air, are automatically opened.

## Cleaning or Replacing the FPD Mirror Metal Plug

Before cleaning or replacing the FPD mirror metal plug, read the following precautions:



**WARNING** Carry out all the operations at low temperature to avoid burns. Therefore, before beginning maintenance, cool the detector to room temperature.



**CAUTION** When handling organic solvents you must take precautions to avoid health hazards.

### Materials needed

FPD fixing tools

Paper towels

Methylene chloride or GC-grade methanol

1-mm Allen wrench


### ❖ To clean or replace the FPD mirror metal plug

1. Put the GC in standby condition.
2. Switch off the flame. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.
6. Remove the FPD detector from the base.
  - a. Disconnect the signal, excitation voltage, and ignition/heating cables from the detector. See [Figure 280](#).

**Figure 280.** FPD Cables Disconnection

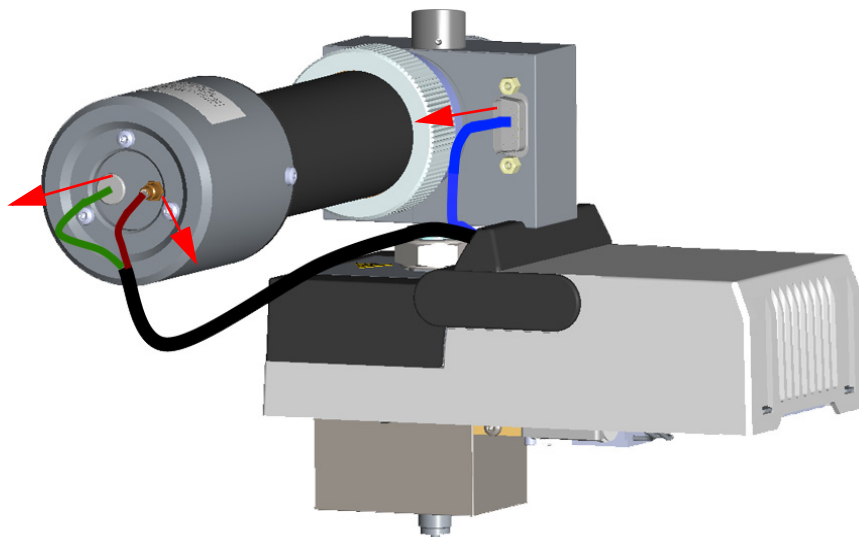

- b. Using the tool provided with the system, loosen the fixing nut on the base of the detector and remove it. See [Figure 281](#).


**Figure 281.** FPD Detector Removal

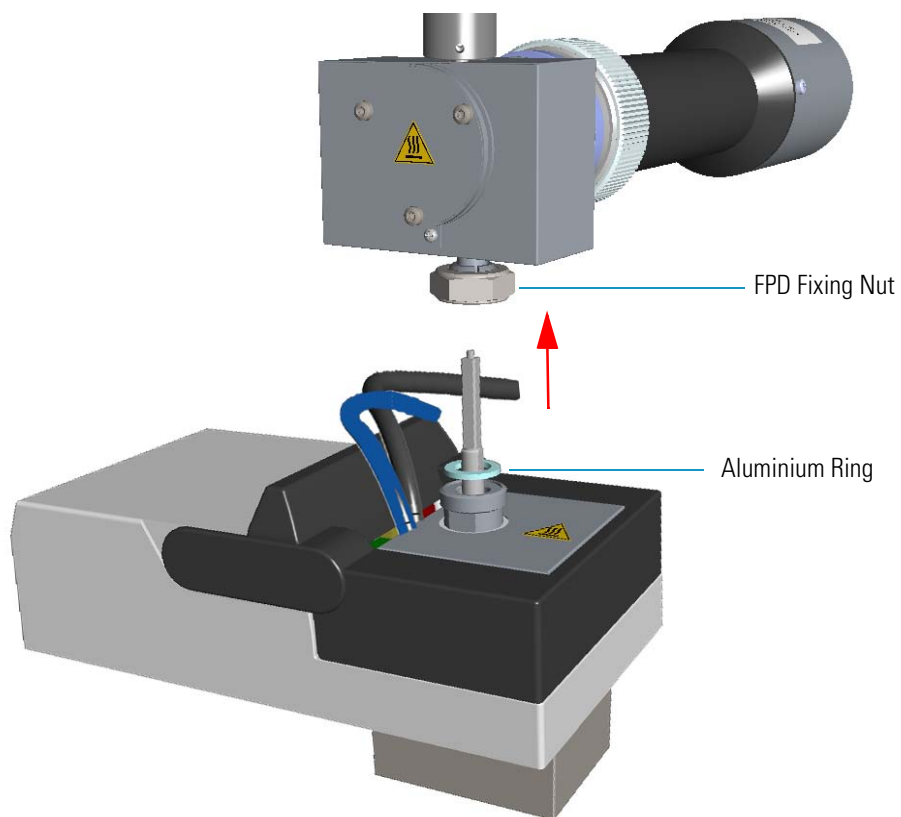

**Note** Do not lose the aluminium ring inserted between the detector head and the base body.

7. Remove the mirror plug.
  - a. Loosen the Allen screws on the mirror plug and remove it. See [Figure 282](#).

**Figure 282.** FPD Mirror Plug Removal

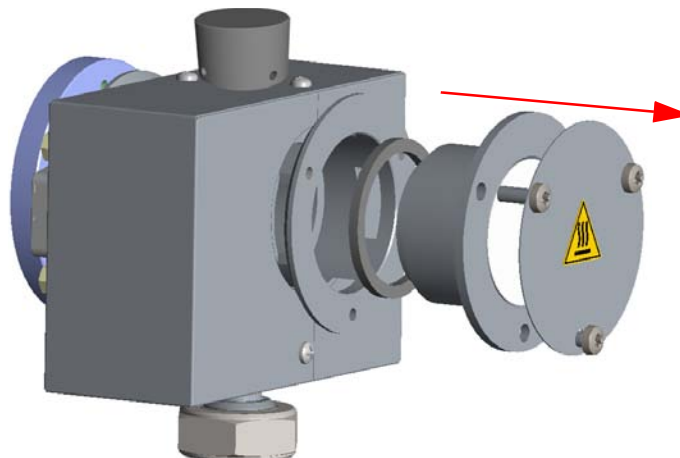


8. Clean the mirror surface of the plug.
  - a. Using a clean paper towel, clean the mirror surface of the plug. If necessary, use a solvent as methylene chloride or methanol to remove deposits, and a metal polishing paste to restore it to the previous reflectivity.
9. Remount the mirror plug.

**Note** Before reinserting the mirror plug, inspect the graphite seal inside the plug's housing. If the seal is damaged and could not ensure tightness, replace it with a new one.

- a. Insert the cleaned or a new mirror plug in the detector body and fix it with the relevant Allen screws. See [Figure 283](#).


**Figure 283.** FPD Mirror Plug Remounting

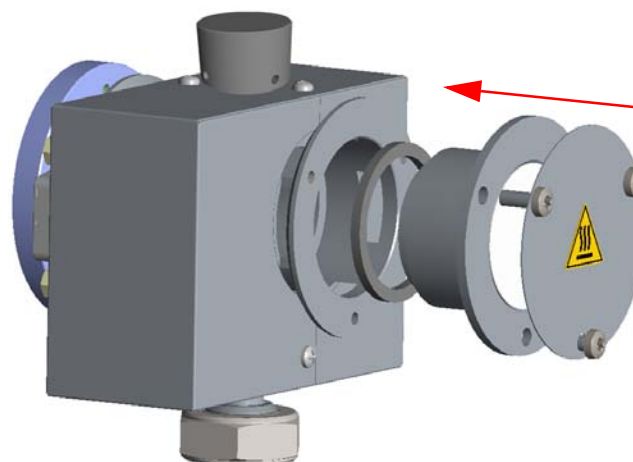

10. Reinstall the FPD detector on the base.


- a. Place the detector on its base body, ensuring that the aluminium ring has been inserted in the correct position, then tighten the fixing nut. See [Figure 284](#).

**Figure 284.** FPD Detector Reinstallation

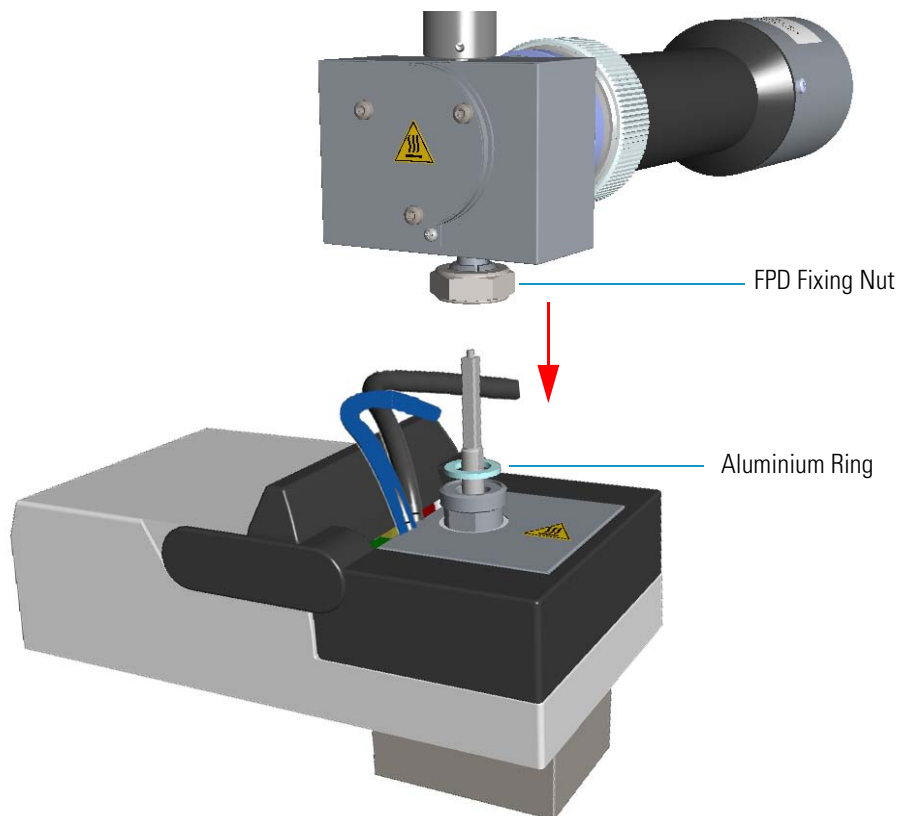

- b. Reconnect the signal, excitation voltage and ignition/heating cables to the detector. See [Figure 285](#).

**Figure 285.** FPD Cables Connection

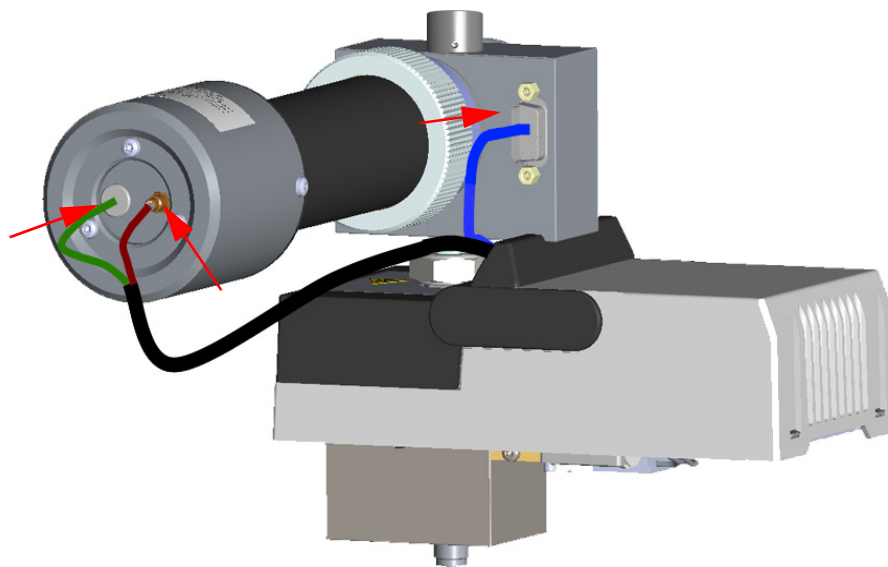



11. If other external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.
12. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
13. Set the normal working conditions.

## Cleaning or Replacing the FPD Filter-side Heat Shields

Before cleaning or replacing the FPD filter-side heat shield, read the following precautions:



**WARNING** Carry out all the operations at low temperature to avoid burns. Therefore, before beginning maintenance, cool the detector to room temperature.



**CAUTION** When handling organic solvents you must take precautions to avoid health hazards.

### Materials needed

FPD fixing tools

Paper towels

FPD maintenance kit

Methylene chloride or GC-grade methanol

1-mm Allen wrench

Screwdriver


### ❖ To clean or replace the FPD filter-side heat shield

1. Put the GC in standby condition.
2. Switch off the flame. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet
5. If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.
6. Remove the FPD detector from the base.
  - a. Disconnect the signal, excitation voltage, and ignition/heating cables from the detector. See [Figure 286](#).

**Figure 286.** FPD Cables Disconnection

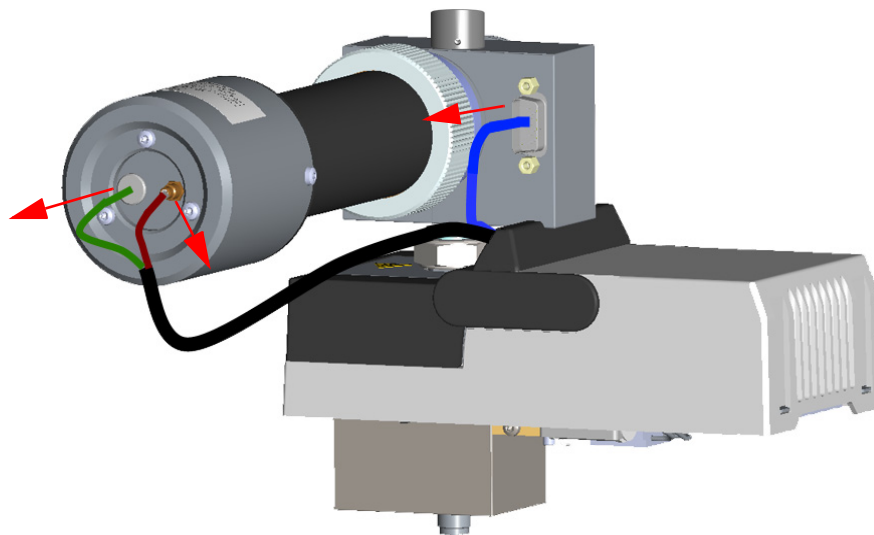

- b. Using the tool provided with the system, loosen the fixing nut on the base of the detector and remove it. See [Figure 287](#).


**Figure 287.** FPD Detector Removal

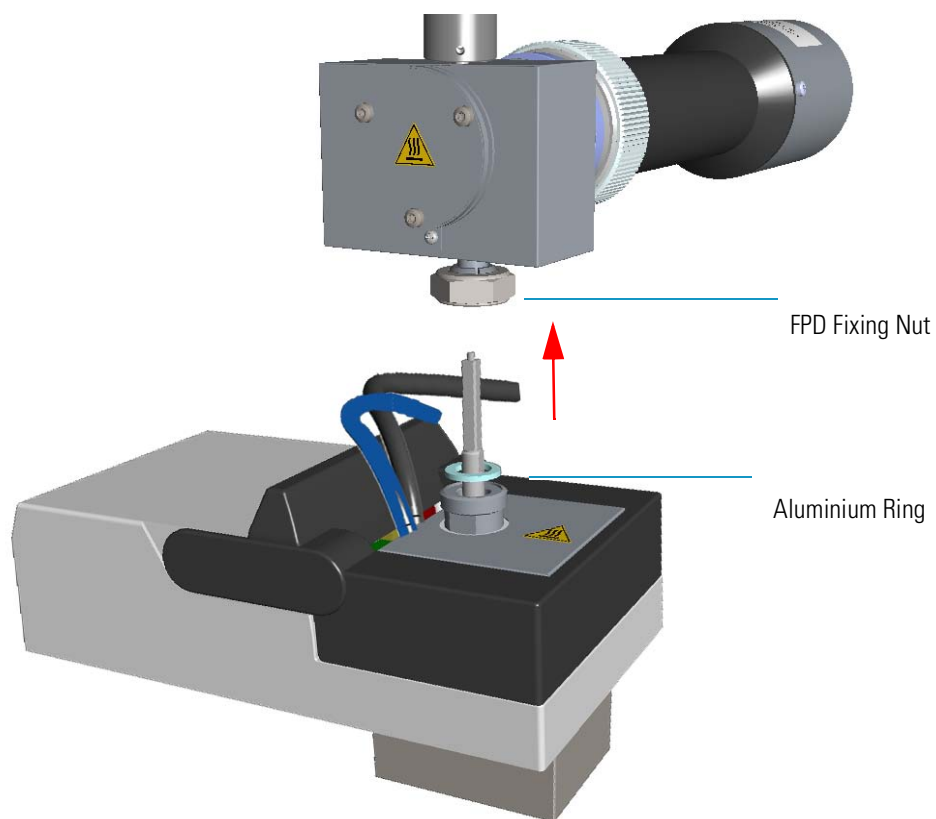

**Note** Do not lose the aluminium ring inserted between the detector head and the base body.

7. Remove the photomultiplier assembly and the filter.
  - a. Loosen the knurled nut that fixes the photomultiplier assembly and remove it from the detector body. See [Figure 288](#).


**Figure 288.** Photomultiplier Assembly Removal

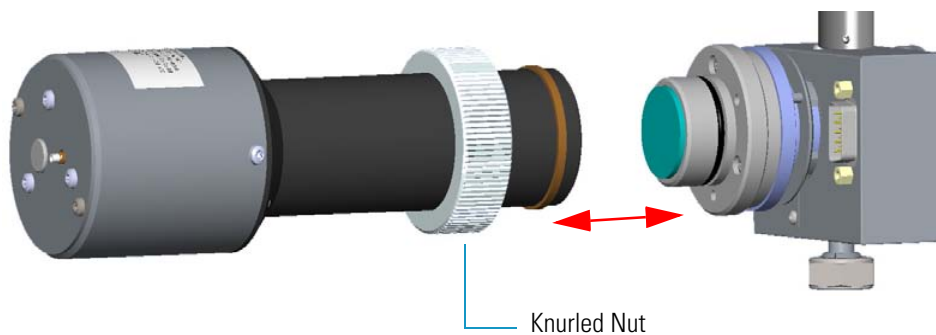

**CAUTION** The photomultiplier tube could damage if exposed to ambient light with the excitation voltage On. Make sure the power supply has been switched off before disconnecting the tube from the detector body.


- b. Remove the interferential filter from its housing, handling it very gently. Keep it using a clean paper towel. See [Figure 289](#).

**Figure 289.** Filter Removal

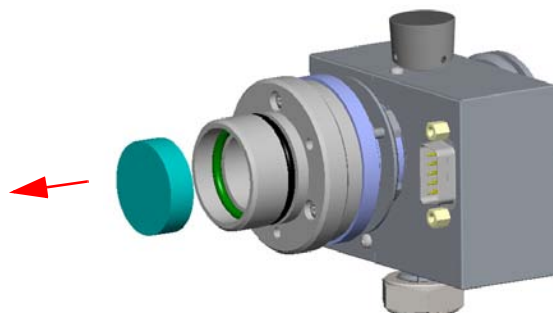

**CAUTION** Filters are fragile. Pay attention not to let the filter fall down and damage.


**CAUTION** If the filter installed is equipped with the spacer (see the kits PN 19050785 and PN 19050786), first remove the spacer, and then the filter.

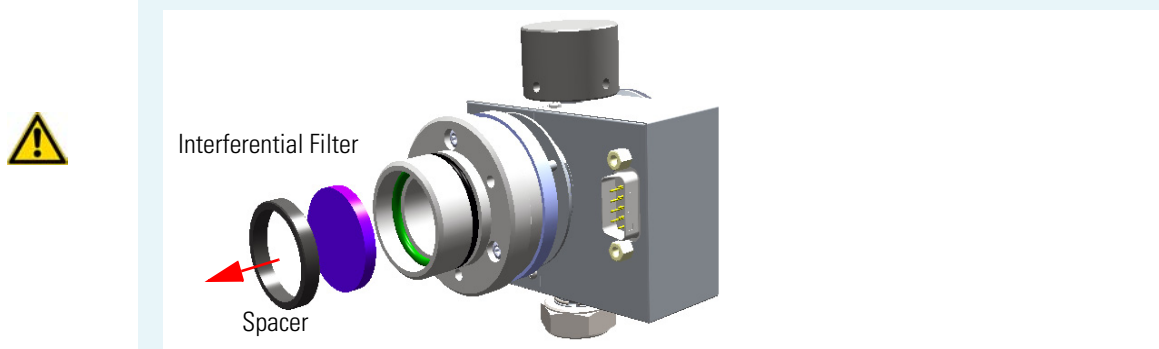

8. Remove the filter-side heat shield.


- a. Loosen the three Allen screws that fix the filter support to the flange of the detector body and remove it. See [Figure 290](#).

**Figure 290.** Filter-side Heat Shield Removal (1)

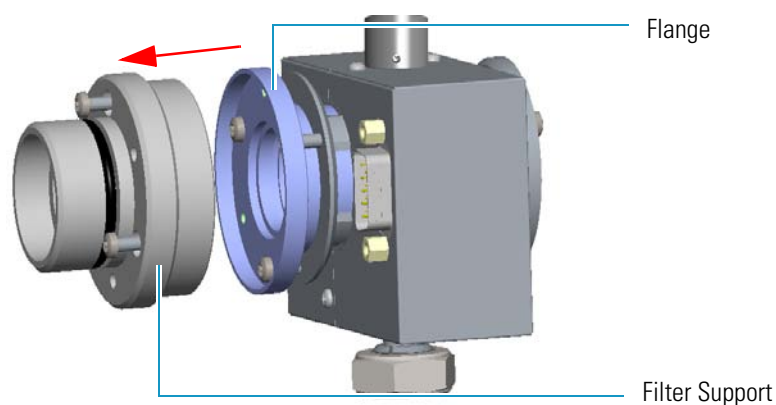

- b. Turn down slowly the detector body and let the filter-side heat shield come out from the flange. Pay attention not to lose the o-ring inserted in the filter support.


**Figure 291.** Filter-side Heat Shield Removal (2)

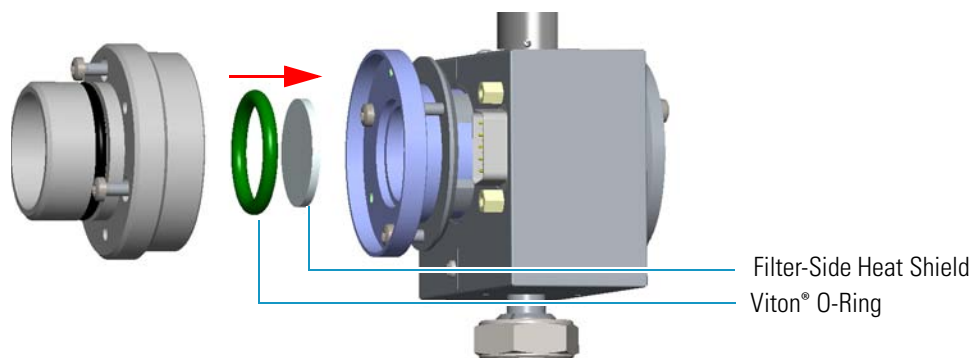

9. Insert a new filter-side heat shield.


- a. Insert a new filter-side heat shield into its housing inside the filter support. See [Figure 292](#).

**Figure 292.** Filter-side Heat Shield Replacement (1)

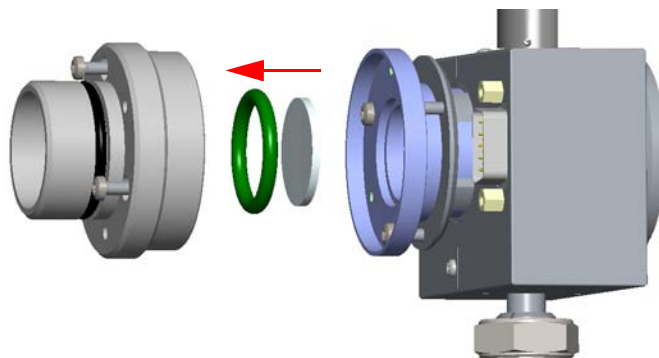

- b. Insert the filter support into the flange and fix it tightening the three Allen screws. Check the o-ring that fixes the heat shield is in place before tightening the screws. See [Figure 293](#).


**Figure 293.** Filter-side Heat Shield Replacement (2)

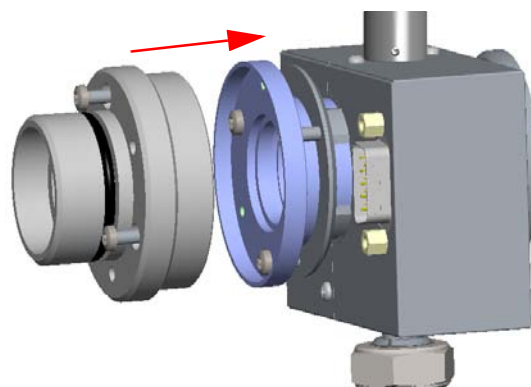

10. Remount the filter and the photomultiplier assembly.

- a. Reinsert the filter into the support. The mirror face must be oriented towards the flame. See [Figure 294](#).


**Figure 294.** Filter Remounting

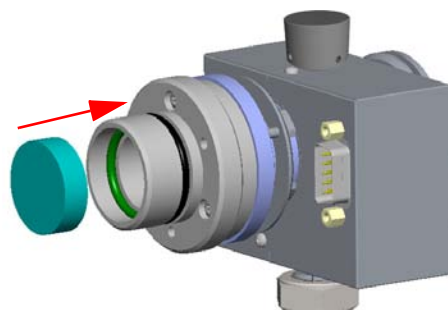




**CAUTION** Avoid touching the filter with your fingers. If you see fingertips on the filter, clean it using a clean paper towel and, if needed, GC-grade methanol before remounting.

**CAUTION** If the filter installed is equipped with the spacer (see the kits PN 19050785 and PN 19050786), first insert the filter, and then the spacer.

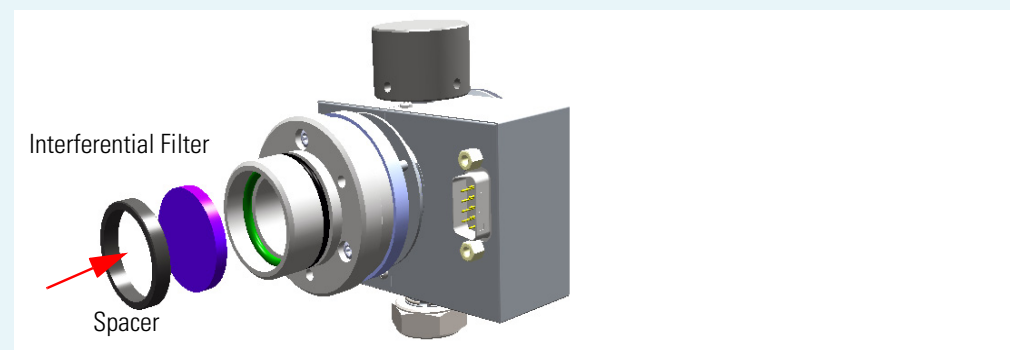


- b. Reassembly the photomultiplier assembly and the detector body, then fix them together tightening the knurled nut. See [Figure 295](#).

**Figure 295.** Photomultiplier Assembly Remounting

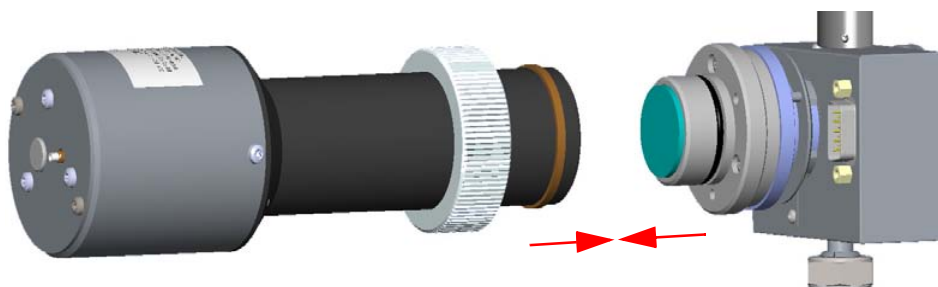


11. Reinstall the FPD detector on the base.
  - a. Place the detector on its base body, ensuring that the aluminium ring has been inserted in the correct position, then tighten the fixing nut. See [Figure 296](#).

**Figure 296.** FPD Detector Reinstallation

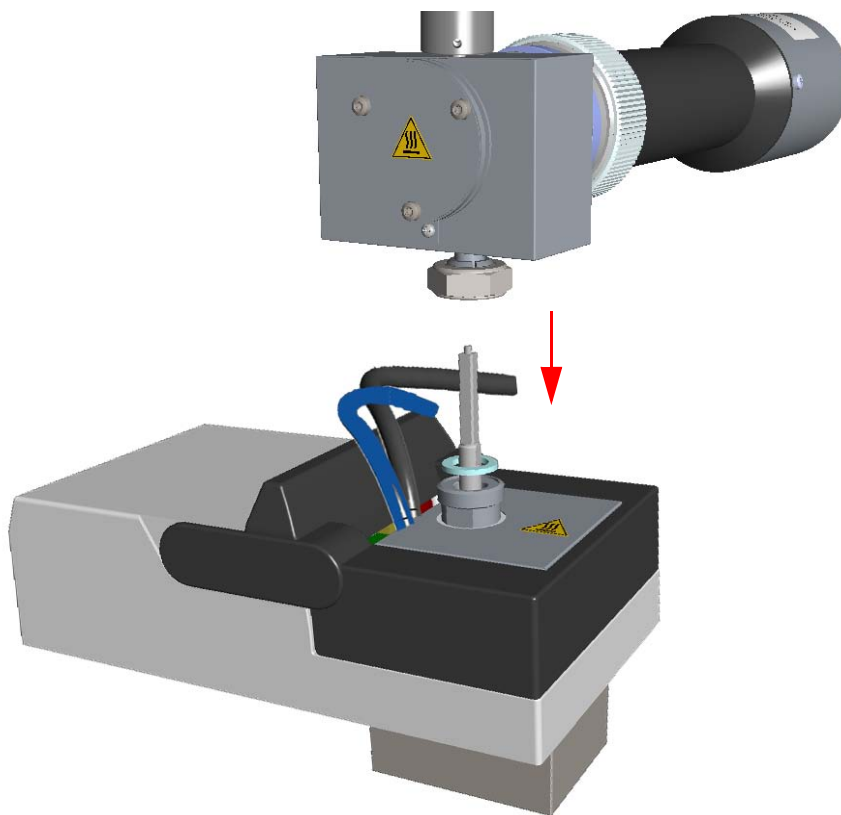

- b. Reconnect the signal, excitation voltage and ignition/heating cables to the detector. See [Figure 297](#).

**Figure 297.** FPD Cables Connection

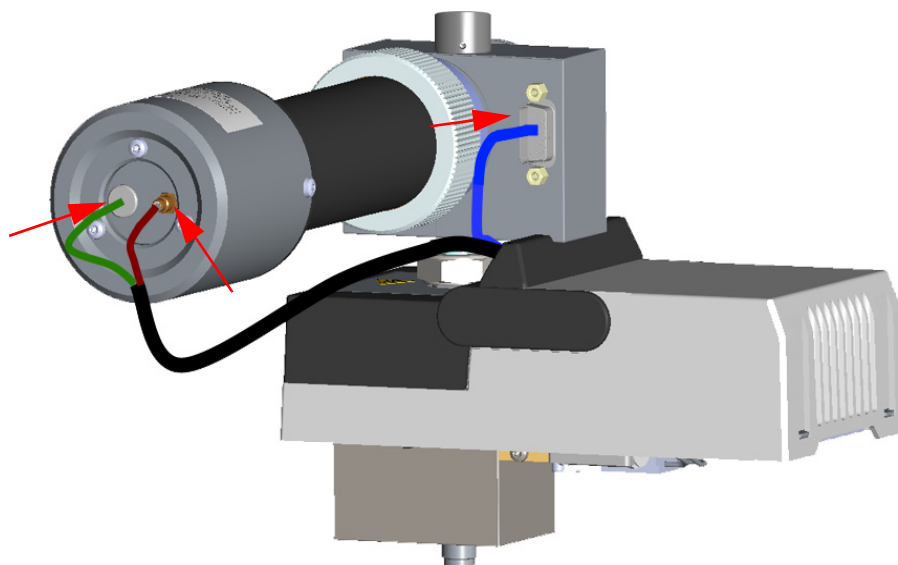

12. If other external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.



13. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
14. Set the normal working conditions.

## Cleaning or Replacing the FPD Flame-side Heat Shields

Before cleaning or replacing the FPD flame-side heat shield, read the following precautions:



**WARNING** Carry out all the operations at low temperature to avoid burns. Therefore, before beginning maintenance, cool the detector to room temperature.



**CAUTION** When handling organic solvents you must take precautions to avoid health hazards.

### Materials needed

FPD fixing tools

Paper towels

FPD maintenance kit

Methylene chloride or GC-grade methanol

1-mm Allen wrench

Screwdriver


### ❖ To clean or replace the FPD flame-side heat shield

1. Put the GC in standby condition.
2. Switch off the flame. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.
6. Remove the FPD detector from the base.
  - a. Disconnect the signal, excitation voltage, and ignition/heating cables from the detector. See [Figure 298](#).

**Figure 298.** FPD Cable Disconnection

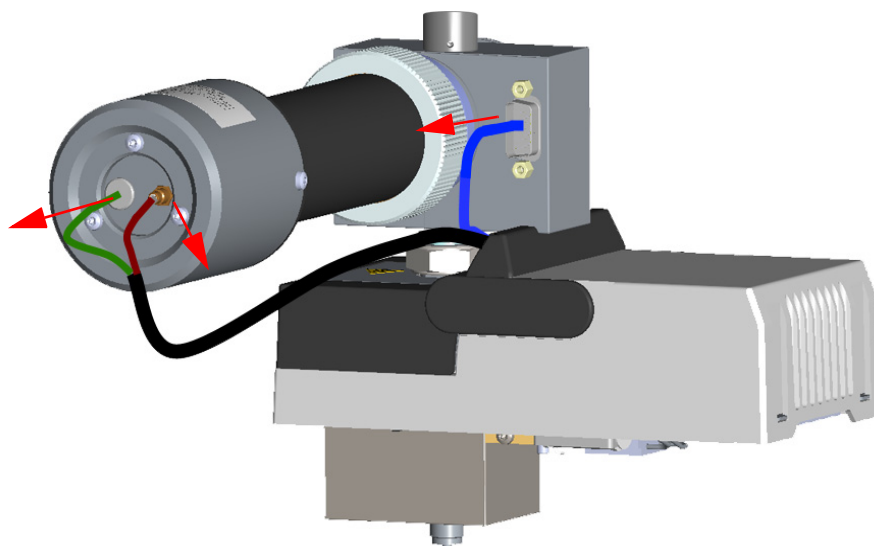

- b. Using the tool provided with the system, loosen the fixing nut on the base of the detector and remove it. See [Figure 299](#).


**Figure 299.** FPD Detector Removal

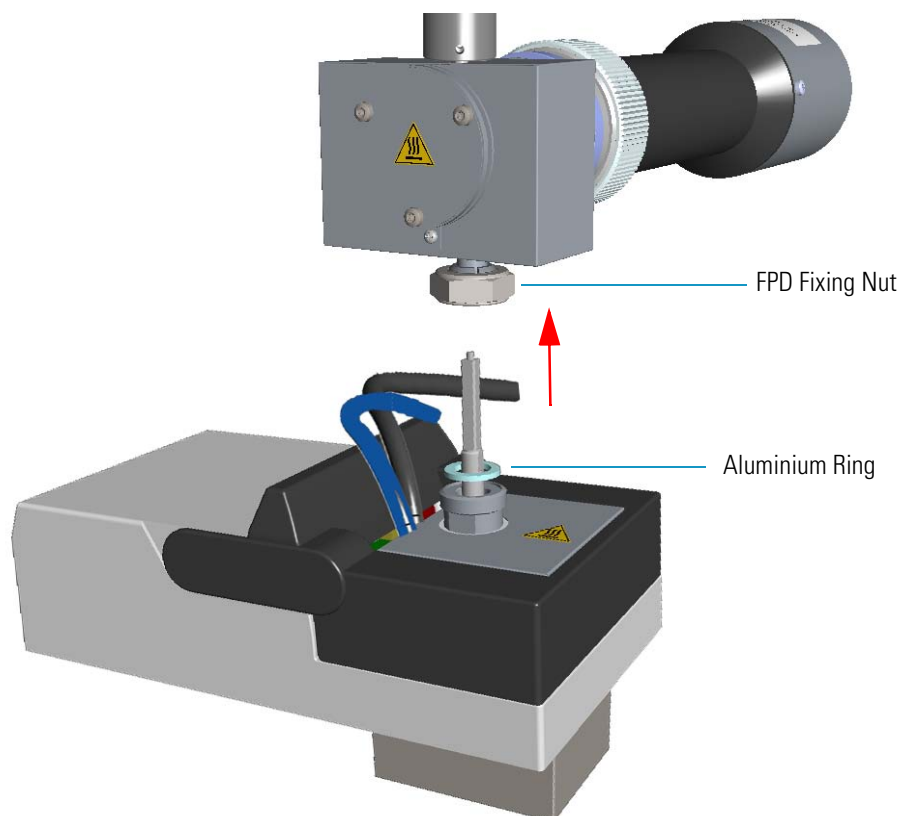

**Note** Do not lose the aluminium ring inserted between the detector head and the base body.

7. Remove the photomultiplier assembly and the filter. See [Figure 300](#).

**Figure 300.** Photomultiplier Assembly Removal

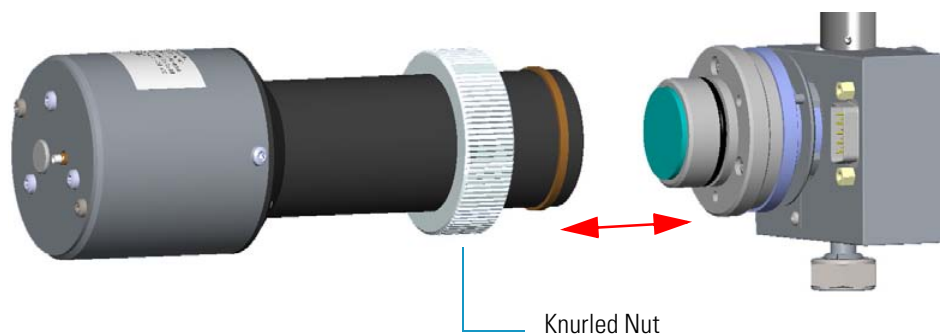


- a. Loosen the knurled nut that fixes the photomultiplier assembly and remove it from the detector body.



**CAUTION** The photomultiplier tube could damage if exposed to ambient light with the excitation voltage On. Make sure the power supply has been switched off before disconnecting the tube from the detector body.


- b. Remove the interferential filter from its housing, handling it very gently. Keep it using a clean paper towel. See [Figure 301](#).

**Figure 301.** Filter Removal

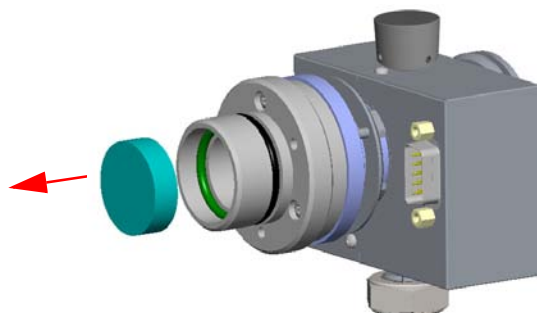

**CAUTION** Filters are fragile. Pay attention not to let the filter fall down and damage.


**CAUTION** If the filter installed is equipped with the spacer (see the kits PN 19050785 and PN 19050786), first remove the spacer, and then the filter.

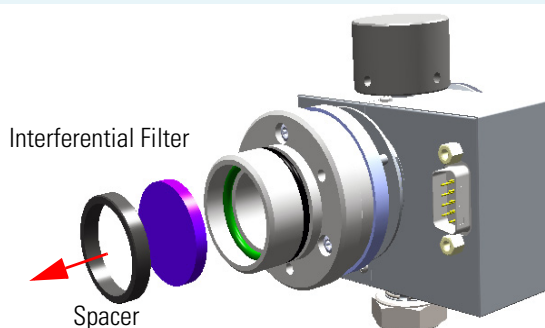

8. Remove the filter support.

- a. Loosen the three Allen screws that fix the filter support to the flange of the detector body and remove it. See [Figure 302](#).


**Figure 302.** Filter Support Removal

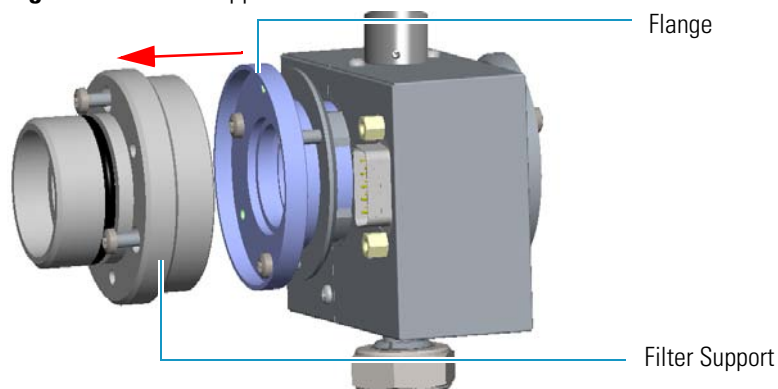

9. Remove the filter-side heat shield.


- a. Loosen the three Allen screws that fix the flange to the detector body and remove it with the spacer that should remain inserted in the flange. Now you could access the flame-side heat shield, and the relevant graphite seal. See [Figure 303](#).

**Figure 303.** Filter-side Heat Shield Removal (1)

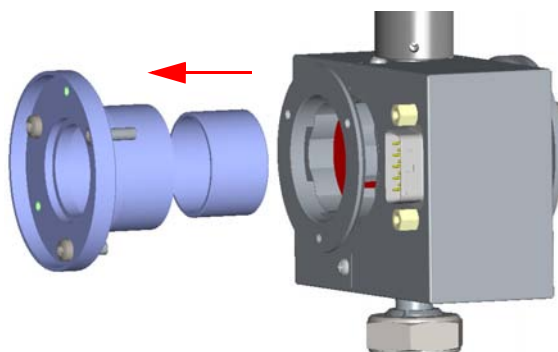

- b. Loosen the Allen screws that fix the mirror plug and remove it. See [Figure 304](#).

**Figure 304.** Mirror Plug Removal

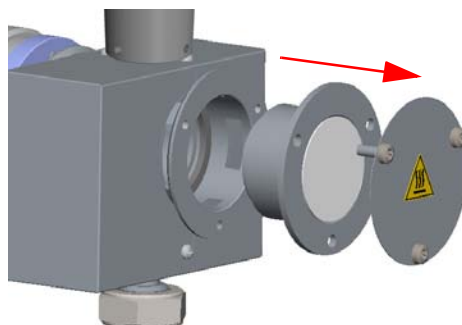


- c. Insert the handle of a screwdriver or other un-sharpened tool in the combustion chamber and push the flame-side heat shield and its graphite seal out from its housing. Act gently to avoid breaking of the heat shield. See [Figure 305](#).




**CAUTION** While pushing out the heat shield, pay attention not to damage the ignition coil.

**Figure 305.** Filter-side Heat Shield Removal (2)

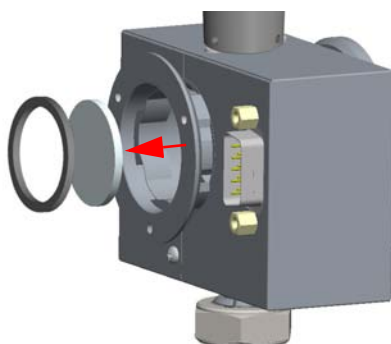


- d. Remove with care traces of graphite due to the breaking of the seal.
10. Insert a new filter-side heat shield. See [Figure 306](#).

**Figure 306.** Filter-side Heat Shield Replacement (1)

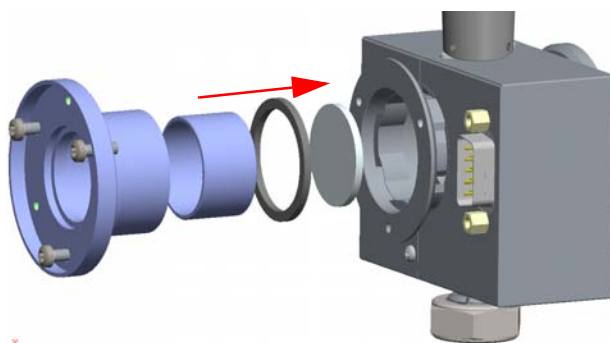

- a. Insert a new heat shield into its housing inside the detector body.
- b. Insert the spacer into the flange, letting it stand out for about 5 mm.
- c. Insert a new graphite seal on the spacer, pushing it slowly until it touches the flange.
- d. Insert the flange into the detector body and fix it tightening the three Allen screws.  
See [Figure 307](#).


**Figure 307.** Filter-side Heat Shield Replacement

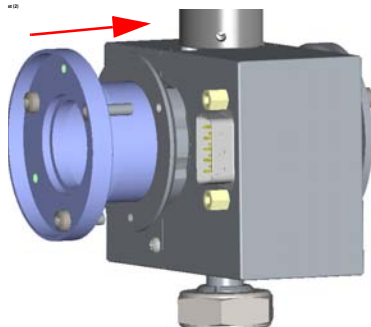

- e. Using a clean paper towel, clean the mirror surface of the plug. If necessary, use a solvent as methylene chloride or methanol to remove deposits and a metal polishing paste to restore it to the previous reflectivity.
- f. Insert the mirror plug in the detector body and fix it with the relevant Allen screws.  
See [Figure 308](#).


**Note** Before reinserting the plug, inspect the graphite seal inside the plug's housing. If the seal is damaged and could not ensure tightness, remove it and replace with a new one.

**Figure 308.** Mirror Plug Reinstallation

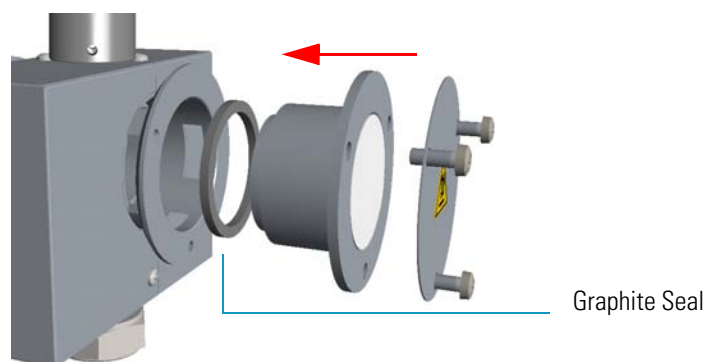

- g. Insert the filter support into the flange and fix it tightening the three Allen screws. Check the O-ring that fixes the heat shield is in place before tightening the screws.  
See [Figure 309](#).


**Figure 309.** Filter Support Reinstallation

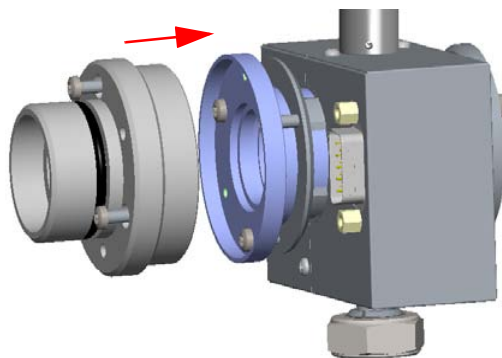

11. Reassembly the filter and the photomultiplier assembly.


- a. Reinsert the filter into the support. The mirror face must be oriented towards the flame. See [Figure 310](#).

**Figure 310.** Filter Remounting

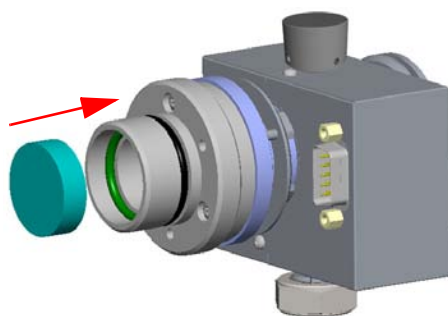

**CAUTION** Avoid touching the filter with your fingers. If you see fingertips on the filter, clean it using a clean paper towel and, if needed, GC-grade methanol before remounting.


**CAUTION** If the filter installed is equipped with the spacer (see the kits PN 19050785 and PN 19050786), first insert the filter, and then the spacer.

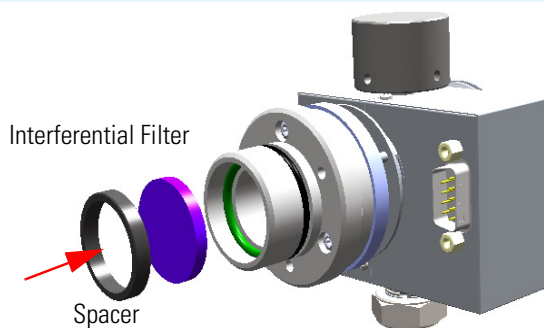

- b. Reassembly the photomultiplier assembly and the detector body, then fix them together tightening the knurled nut. See [Figure 311](#).




**Figure 311.** Photomultiplier Assembly Remounting

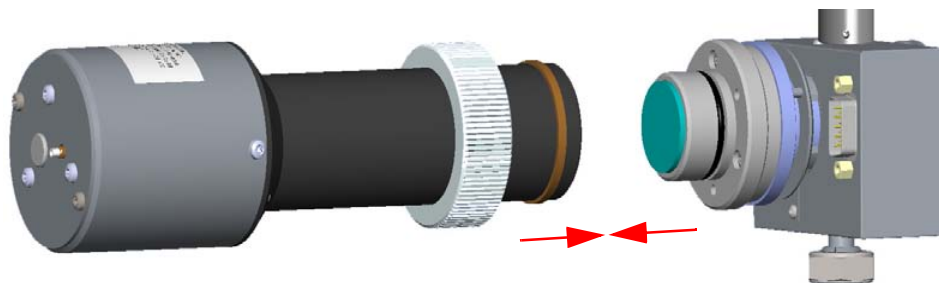

12. Reinstall the FPD detector on the base.


- a. Place the detector on its base body, ensuring that the aluminium ring has been inserted in the correct position, then tighten the fixing nut. See [Figure 312](#).

**Figure 312.** FPD Detector Reinstallation

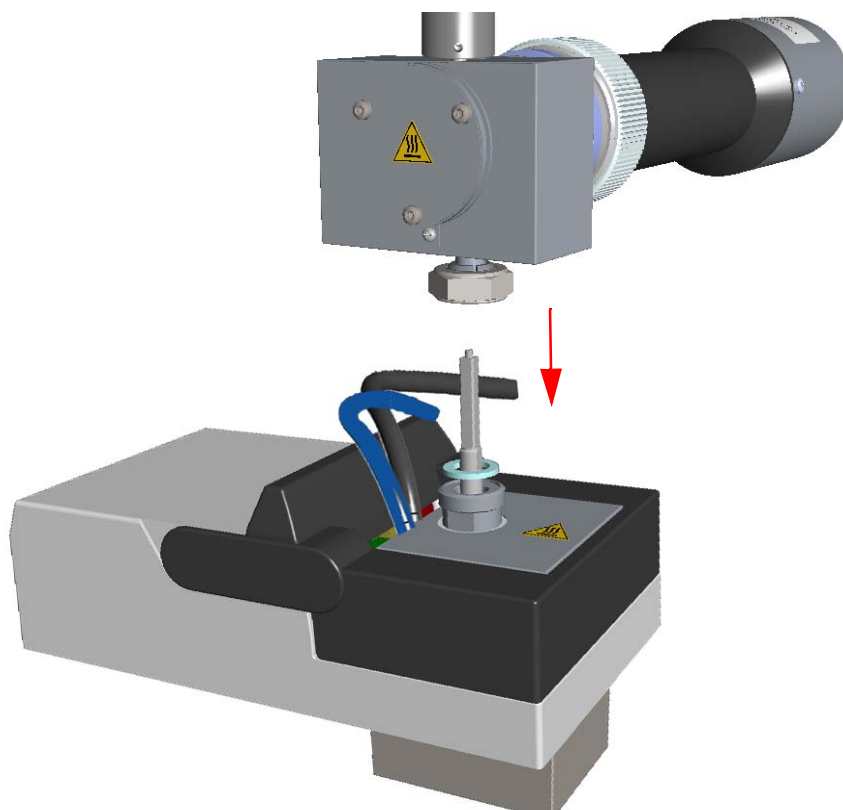

- b. Reconnect the signal, excitation voltage and ignition/heating cables to the detector. See [Figure 313](#).

**Figure 313.** FPD Cables Connection

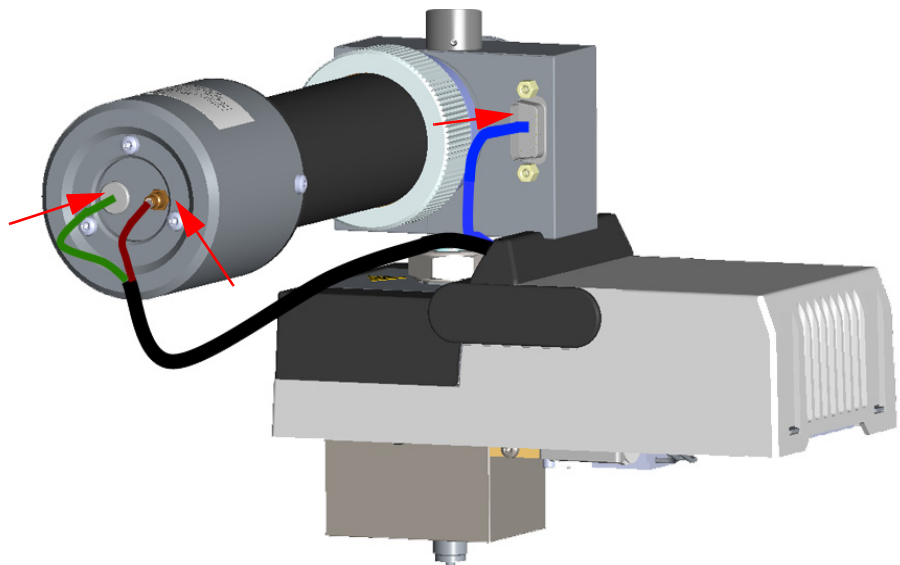

13. If other external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.
14. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
15. Set the normal working conditions.

# Replacing the FPD Photomultiplier Tube

Before replacing the FPD photomultiplier tube, read the following precautions:



**WARNING** Carry out all the operations at low temperature to avoid burns. Therefore, before beginning maintenance, cool the detector to room temperature.



**CAUTION** The photomultiplier tube must be replaced only when defective. Wear clean, lint- and powder-free gloves when you handle the photomultiplier tube.

## Materials needed

FPD fixing tools

FPD Photomultiplier Tube

1-mm Allen wrench

Cross head screwdriver

Gloves


## ❖ To replace the photomultiplier tube

1. Put the GC in standby condition.
2. Switch off the flame. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.
6. Remove the FPD detector from the base.
  - a. Disconnect the signal, excitation voltage, and ignition/heating cables from the detector. See [Figure 314](#).

**Figure 314.** FPD Cable Disconnection

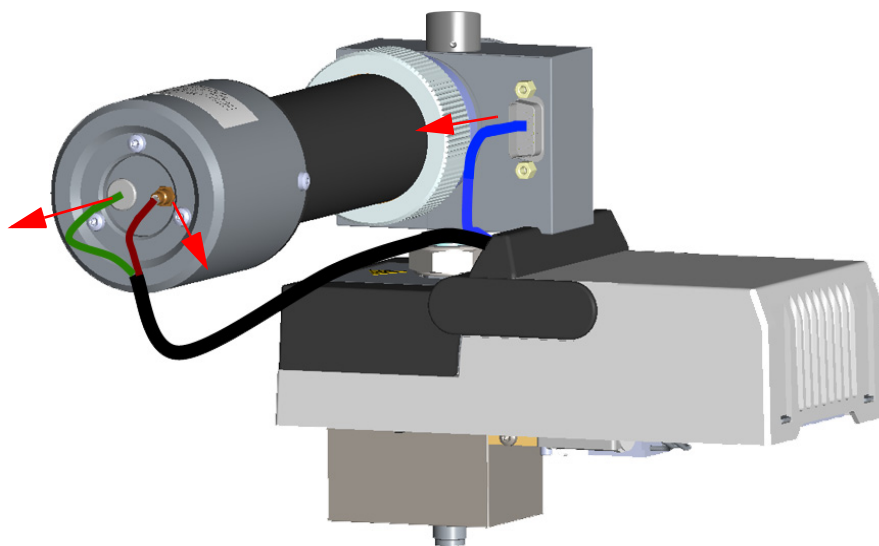

- b. Using the tool provided with the system, loosen the fixing nut on the base of the detector and remove it. See [Figure 315](#).


**Figure 315.** FPD Detector Removal

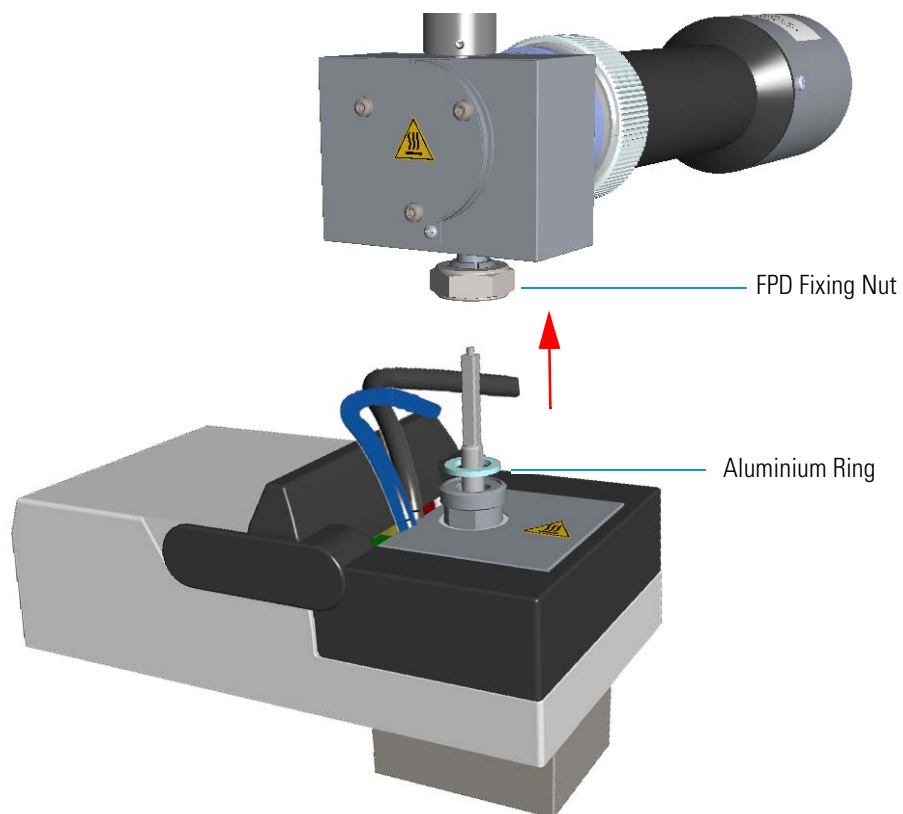

**Note** Do not lose the aluminium ring inserted between the detector head and the base body.


7. Dismount the photomultiplier tube assembly.
  - a. Using a crosshead screwdriver unscrew the two fixing screws that fix the photomultiplier tube housing to the photomultiplier tube holder. See [Figure 316](#).

**Figure 316.** Photomultiplier Tube Housing Removal (1)

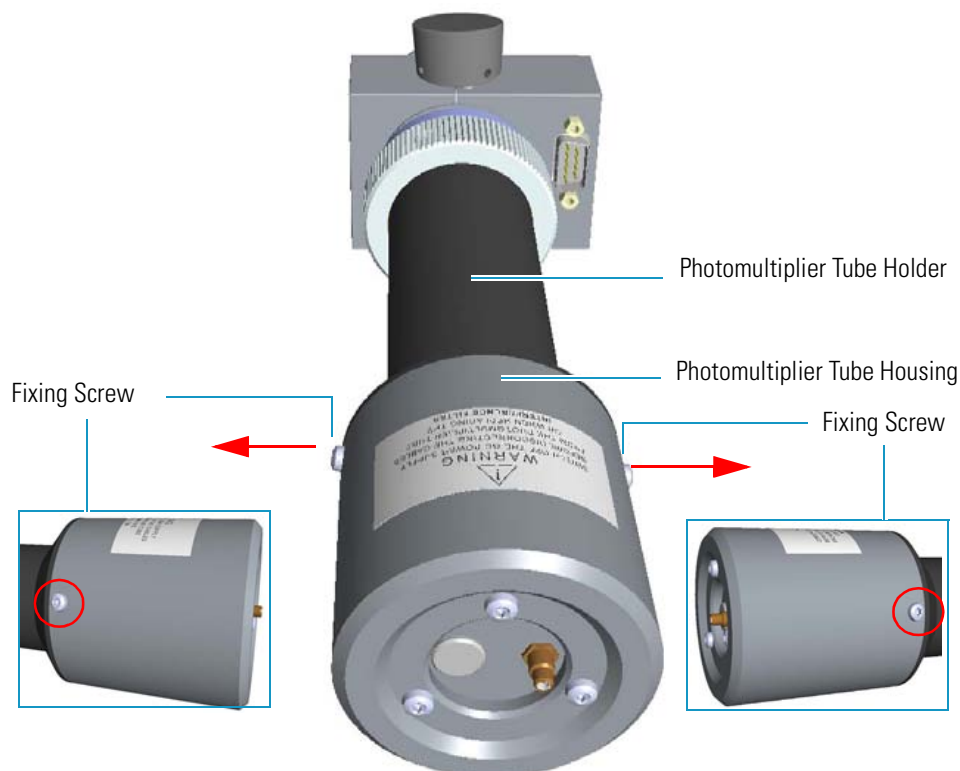

- b. Remove the photomultiplier tube housing from the photomultiplier tube holder. Holding the photomultiplier tube holder with one hand, use the other hand to pull out the photomultiplier tube housing from the holder. During the extraction meets some resistance due to the presence of an o-ring inside the holder. See [Figure 317](#).

**Figure 317.** Photomultiplier Tube Housing Removal (2)

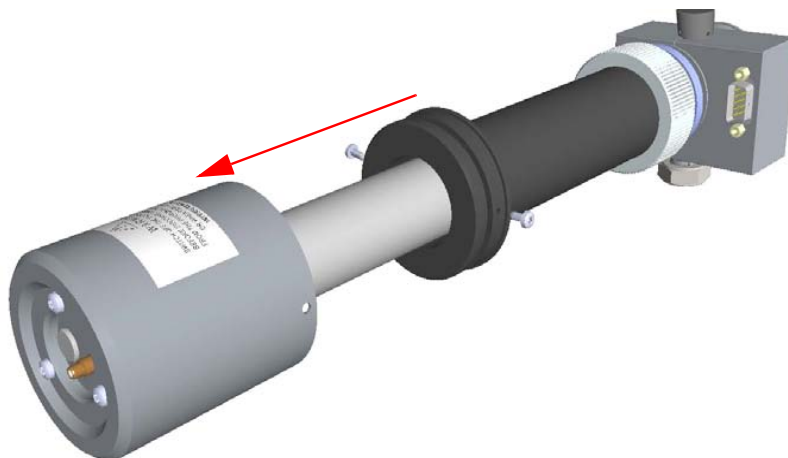


8. Replace the defective photomultiplier tube.

**Note** Wear clean, lint- and powder-free gloves when you handle the photomultiplier tube.


- a. Carefully extract the photomultiplier tube from its connector. See [Figure 318](#).

**Figure 318.** Photomultiplier Tube Replacement (1)

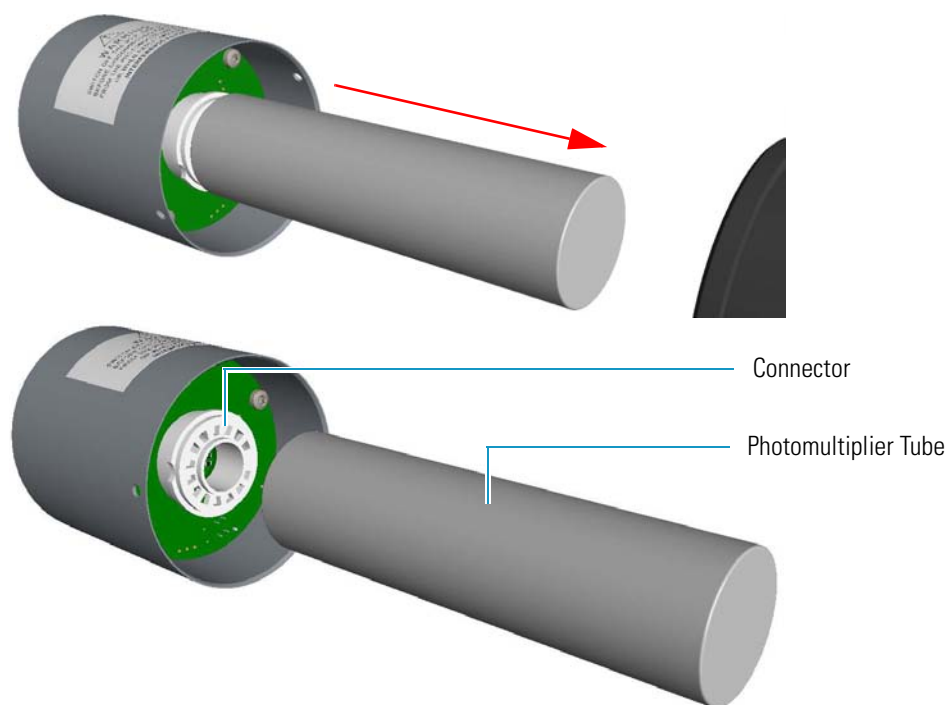


- b. Replace the photomultiplier tube with a new one. Make sure to correctly insert the pins of the photomultiplier tube into the connector referring to the polarized position. See [Figure 319](#).

**Figure 319.** Connector Polarized Position

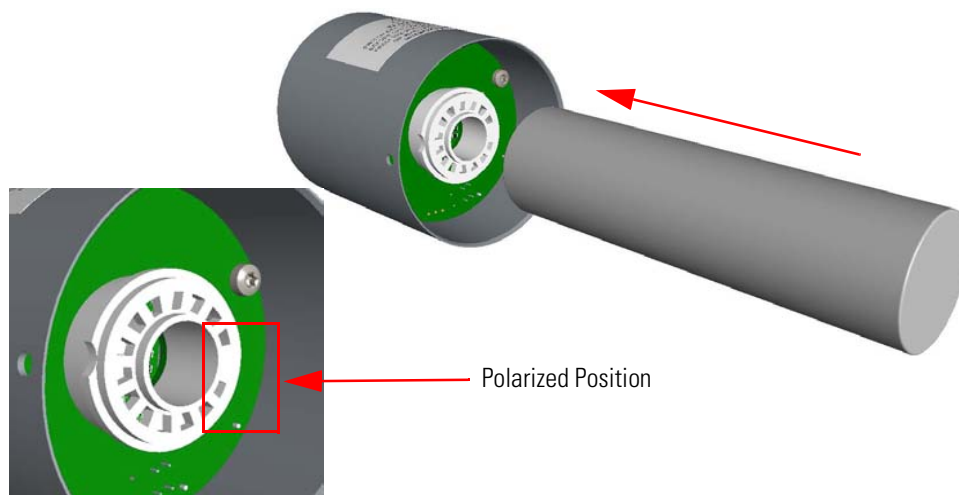


9. Remount the photomultiplier tube assemble.
- Reinsert the photomultiplier tube housing into the photomultiplier tube holder. Holding the photomultiplier tube holder with one hand, use the other hand to reinsert the photomultiplier tube housing into the holder. During the insertion meets some resistance due to the presence of an o-ring inside the holder. See [Figure 320](#).

**Figure 320.** Photomultiplier Tube Housing Reinstallation (1)

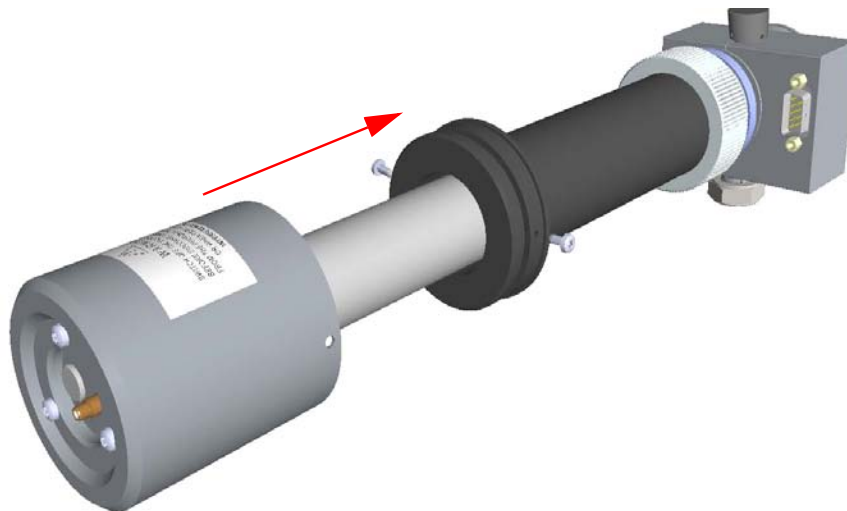


- Fix the photomultiplier tube housing to the photomultiplier tube holder. Carefully rotate the photomultiplier tube housing up to matching its fixing holes to the fixing holes on the photomultiplier tube holder. See [Figure 321](#).

**Figure 321.** Photomultiplier Tube Housing Reinstallation (2)

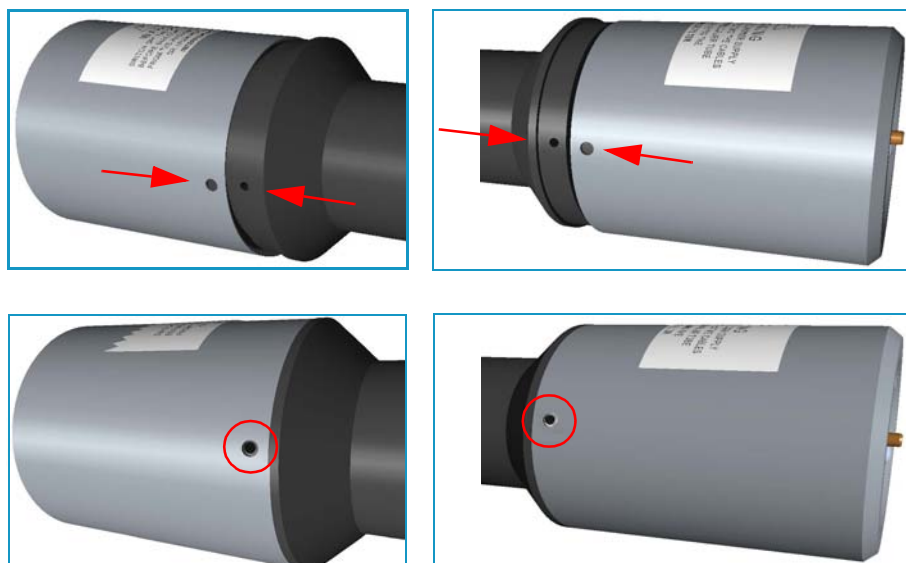

- Insert the two fixing screws into the corresponding fixing holes, then tighten the screws using a crosshead screwdriver.


**Figure 322.** Photomultiplier Tube Housing Reinstallation (3)

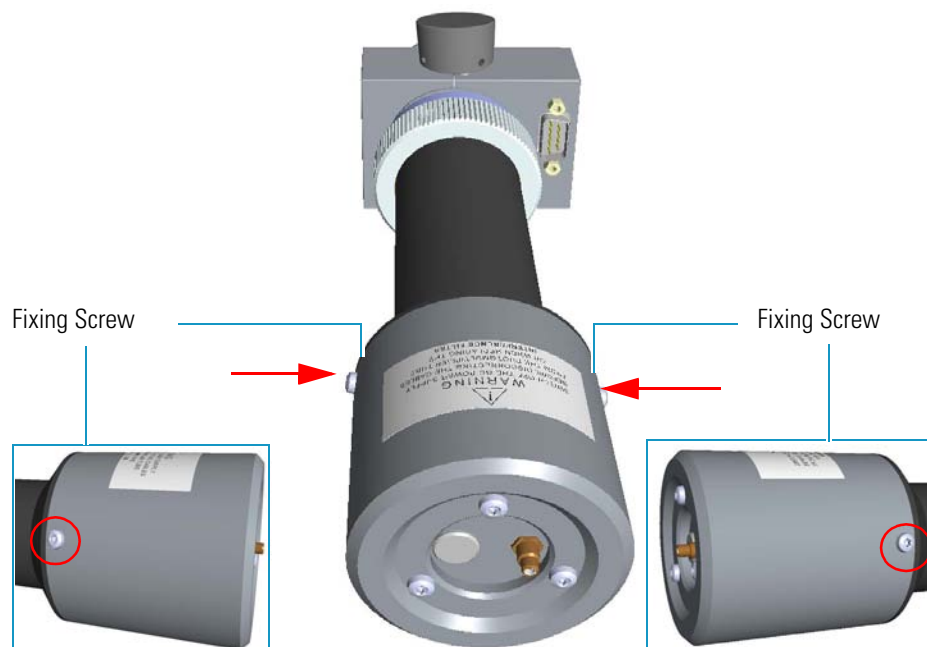

10. Reinstall the FPD detector on the base.

- a. Place the detector on its base body, ensuring that the aluminium ring has been inserted in the correct position, then tighten the fixing nut. See [Figure 323](#).


**Figure 323.** FPD Detector Reinstallation

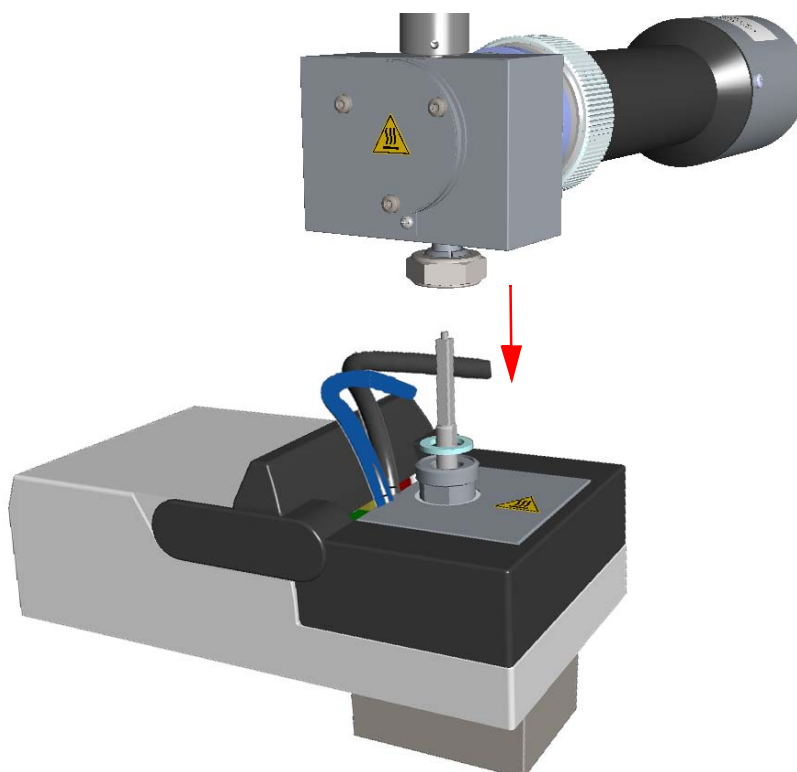



- b. Reconnect the signal, excitation voltage and ignition/heating cables to the detector.  
See [Figure 324](#).

**Figure 324.** FPD Cables Connection

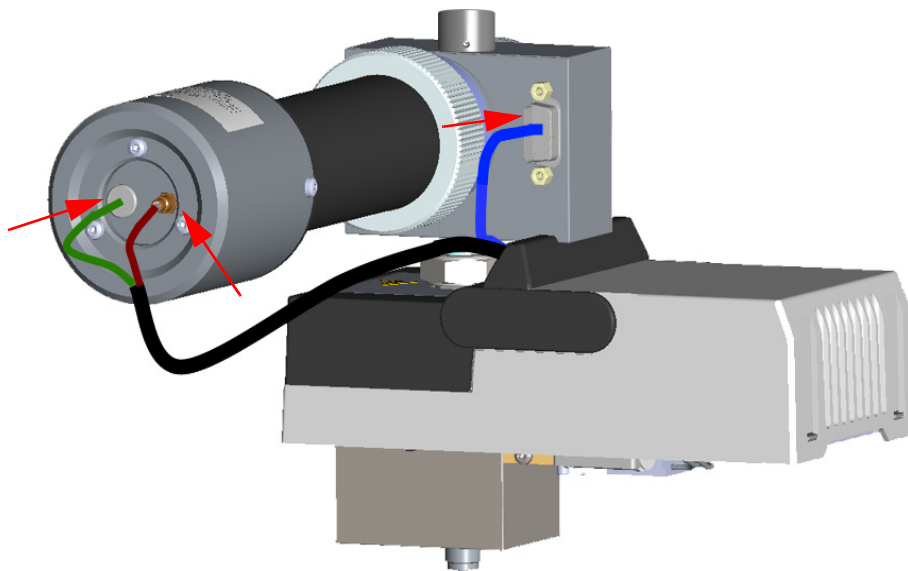

11. If other external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.
12. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
13. Set the normal working conditions.

## **7 Detectors Advanced Maintenance**

Replacing the FPD Photomultiplier Tube

## Installing Optional Kits

This chapter describes how to install the optional kits available for the TRACE 1300/TRACE 1310. See the *TRACE 1300 and TRACE 1310 Spare Parts Guide* for information about ordering the upgrade kits in this chapter.

### Contents

- Installing the Oven Exhaust Kit
- Installing the Merlin Microseal High Pressure Valve Kit
- Installing the Purge & Trap Adapter Kit on the SSL/SSLBKF Injector
- Installing the Packed Column Adapters
- Installing the HS Adapter Kit on the SSL/SSLBKF Injector
- Installing the Large Volume Splitless Kit
- Installing the Manual On/Off Valve for Single Gas Line
- Connecting a SSL/PTV Backflush System
- Connecting a SSL/PTV Backflush System for High Temperature
- Connecting a GSV Backflush System
- Installing the NoVent Microfluidics
- Installing a FTIR Make-up Module
- Installing the Hot Injection Adapter Kit on the SSL/SSLBKF Injector
- Performing the Dual FPD Detector Configuration

## Installing the Oven Exhaust Kit

The oven vents at the back of the GC discharge hot air up to 450 °C (842 °F) during cooling. Installing the optional Oven Exhaust Kit (PN 19050760), the hot air from the oven vents can be carried away to a fume hood or other exhausting devices.


The kit includes two air ducts and two 3-m length extensible aluminium tubes. See [Figure 325](#).

**Figure 325.** Air Duct and Extensible Aluminium Tube



❖ **To install the oven exhaust kit**

1. Put the GC in standby condition.
2. Cool the oven to room temperature.
3. Place each air duct on the oven vent inserting the hooks of the air duct into the relevant vertical slots provided on the GC back panel.
4. Place an end of each extensible aluminium tube on each air duct.
5. Connect the other end of the extensible aluminium tube to an exhausting device.

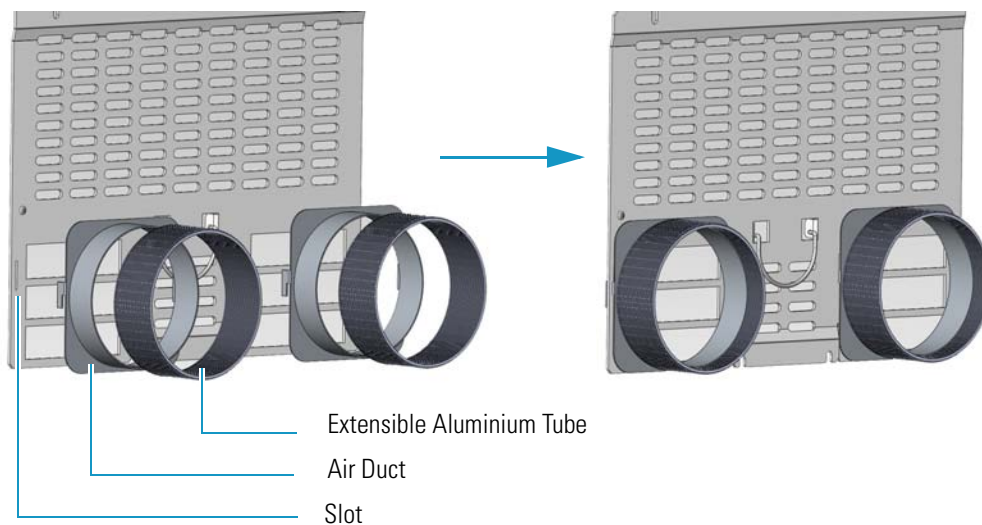

6. Set the normal GC working conditions.

## Installing the Merlin Microseal High Pressure Valve Kit

This section provides instructions for installing the Merlin Microseal™ High Pressure Valve kit (PN 19050735) on the SSL, SSLBKF, HeS-S/SL, PTV, and PTVBKF injectors.

**Note** Merlin Microseal™ is a trademark of the Merlin Instrument Company.

### Introduction

The Merlin Microseal™ High Pressure Valve is a long-life replacement for the standard septum on the SSL, SSLBKF, HeS-S/SL, PTV, and PTVBKF injectors of the TRACE 1300/TRACE 1310. High pressure capability allows operation from 15 to 700 kPa (2 - 100 psi). A top wiper rib improves resistance to particulate contamination. The valve can be taken apart for cleaning. High resistance to wear greatly reduces the shedding of valve particles into the injection port liner. This eliminates a major source of septum bleed and ghost peaks in the chromatogram. Longer life reduces the chances of valve leaks occurring during extended automated runs. The low syringe insertion force makes manual injections easier. The Microseal High Pressure valve requires a 0.63-mm diameter (0.025-in.) blunt tip syringe.

**Syringe Compatibility** — The Microseal valve should be used only with a blunt-tipped 0.63-mm diameter (0.025-in.; 23 gauge) syringe needle.

- A needle with too small a diameter, such as a 0.43-mm diameter (0.017-in.; 23 gauge) needle, will not seal properly.
- A needle with too large a diameter will overstretch and damage the O-ring and duckbill seals.
- A sharp-pointed or sharp-edged syringe needle will slice or pierce the seals.

**Temperature and Pressure Limits** — The operating temperature of the valve is lower than the injection port temperature set point. For long term operation (> 6 months) use the following limits for operating conditions:

- Injection port temperature < 350 °C; Pressure range: 15 - 700 kPa (2 - 100 psi)

Higher temperatures and pressures result in shorter lifetime. High temperature deterioration can be recognized by leaks caused by stiffening and cracking of the Microseal valve, particularly around the sealing flange.

### Getting Started

Install the Microseal High Pressure Valve on the SSL, SSLBKF, HeS-S/SL, PTV, and PTVBKF injector following the instructions in the next operating procedures.

- [“To install the Merlin Microseal Valve kit on the SSL/SSLBKF and HeS-S/SL injector” on page 328](#)
- [“To install the Merlin Microseal Valve kit on the PTV/PTVBKF injector” on page 329](#)

## 8 Installing Optional Kits

### Installing the Merlin Microseal High Pressure Valve Kit


#### ❖ To install the Merlin Microseal Valve kit on the SSL/SSLBKF and HeS-S/SL injector

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


4. Put the autosampler away if present.
5. Open the module flap cover.
6. Remove the septum.

**Figure 326.** SSL/SSLBKF and HeS-S/SL Injector: Septum Replacement

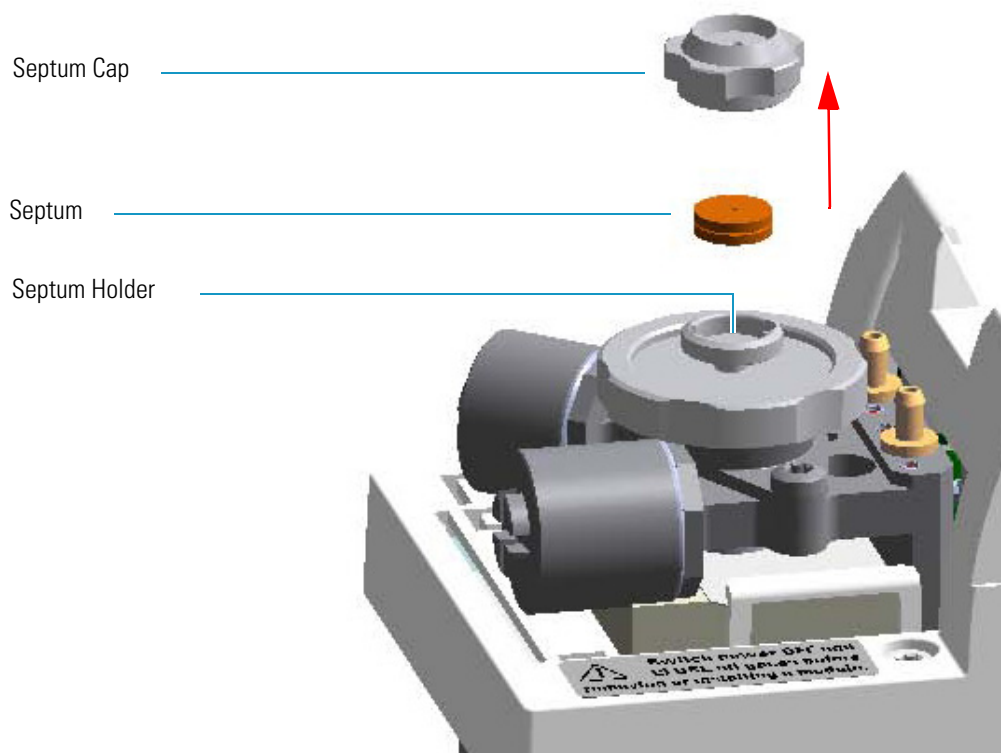

- a. Unscrew and remove the septum cap.
  - b. Using tweezers, remove the septum from the septum holder.
7. Install the Microseal valve on the injector.

**Figure 327.** SSL/SSLBKF and HeS-S/SL Injector: Microseal Valve Installation

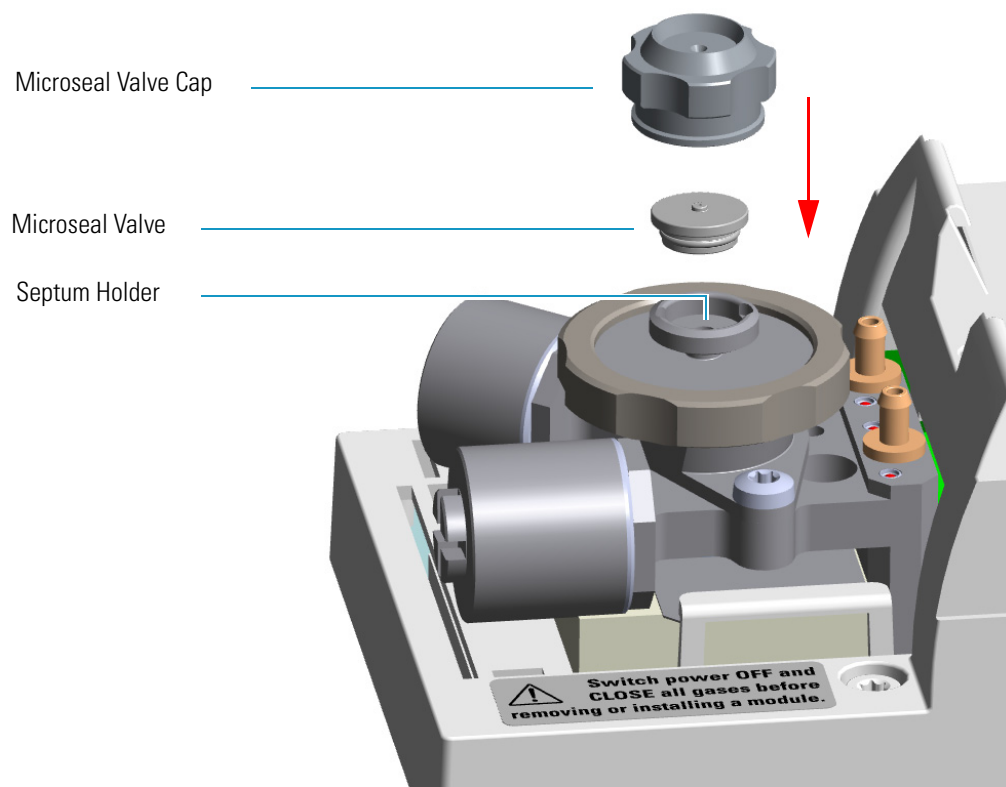

- a. Place the Microseal valve in the septum holder.
  - b. Push down gently to seat the Microseal valve onto the rim of the septum holder.
  - c. Screw the Microseal valve cap completely. Microseal valve overtightening and sealing flange stressing is not mechanically allowed.
8. Close the module flap cover.
  9. If present, move the autosampler towards the module to restore the original alignment.
  10. Turn the carrier gas on.
  11. Set the normal injector, detector, and GC working conditions.

❖ **To install the Merlin Microseal Valve kit on the PTV/PTVBKF injector**

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


4. Put the autosampler away if present.

## 8 Installing Optional Kits

### Installing the Merlin Microseal High Pressure Valve Kit

5. Open the module flap cover.
6. Remove the septum.


**Figure 328.** PTV/PTVBKF Injector: Septum Replacement

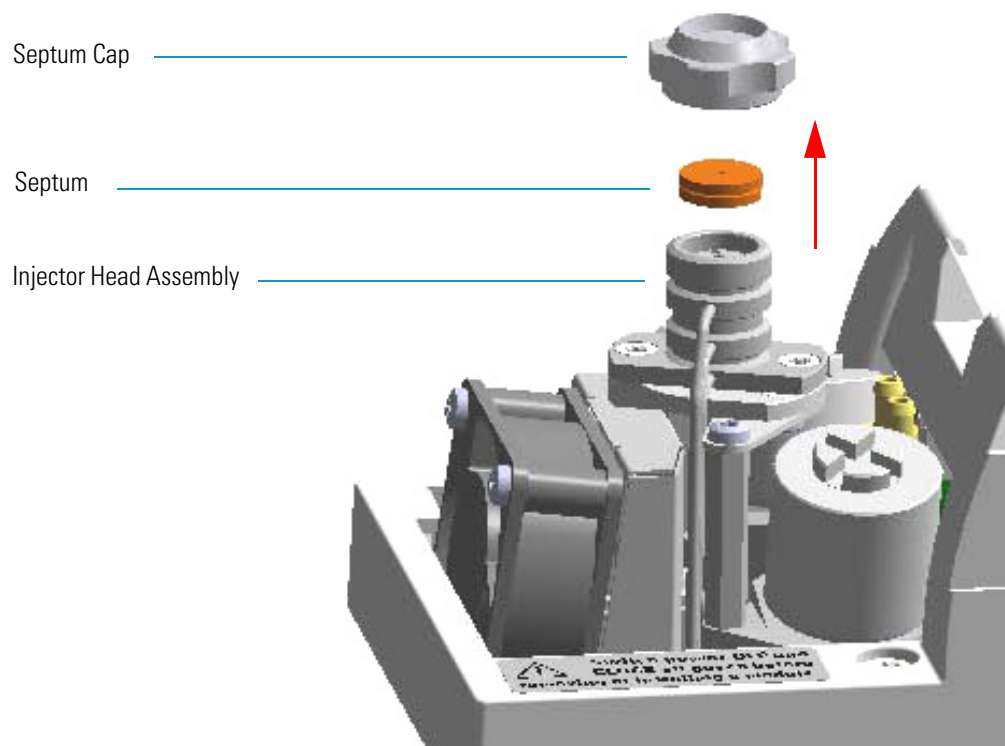

- a. Unscrew the septum cap of the injector.
  - b. Using tweezers, remove the septum from the injector head assembly.
7. Install the Microseal valve on the injector.



**Figure 329.** PTV/PTVBKF Injector: Microseal Valve Installation

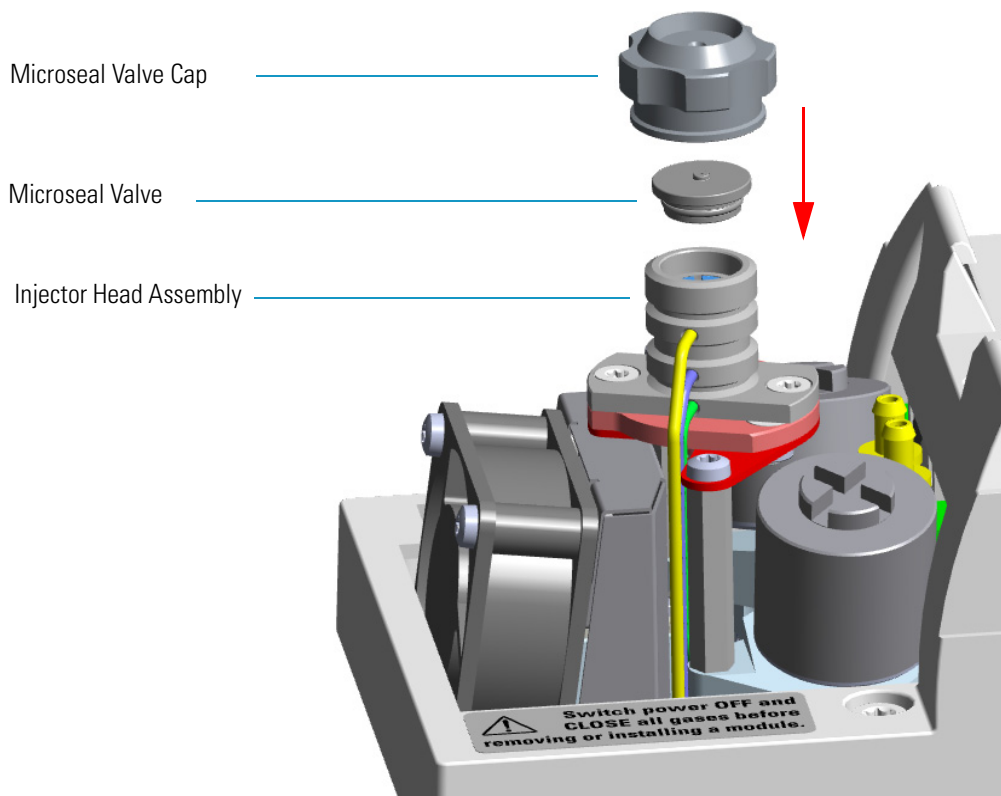


- a. Place the Microseal valve in the injector head assembly.
  - b. Push down gently to seat the Microseal valve onto the rim of the septum holder.
  - c. Screw the Microseal valve cap completely. Microseal valve overtightening and sealing flange stressing is not mechanically allowed.
8. Close the module flap cover.
  9. If present, move the autosampler towards the module to restore the original alignment.
  10. Turn the carrier gas on.
  11. Set the normal injector, detector, and GC working conditions.

# Installing the Purge & Trap Adapter Kit on the SSL/SSLBKF Injector

This section provides the instruction for installing the Purge & Trap Adapter kit (PN 19050730) on the SSL/SSLBKF injector for the use with a Purge & Trap system. See [Figure 330](#).

**Figure 330.** Purge & Trap Adapter

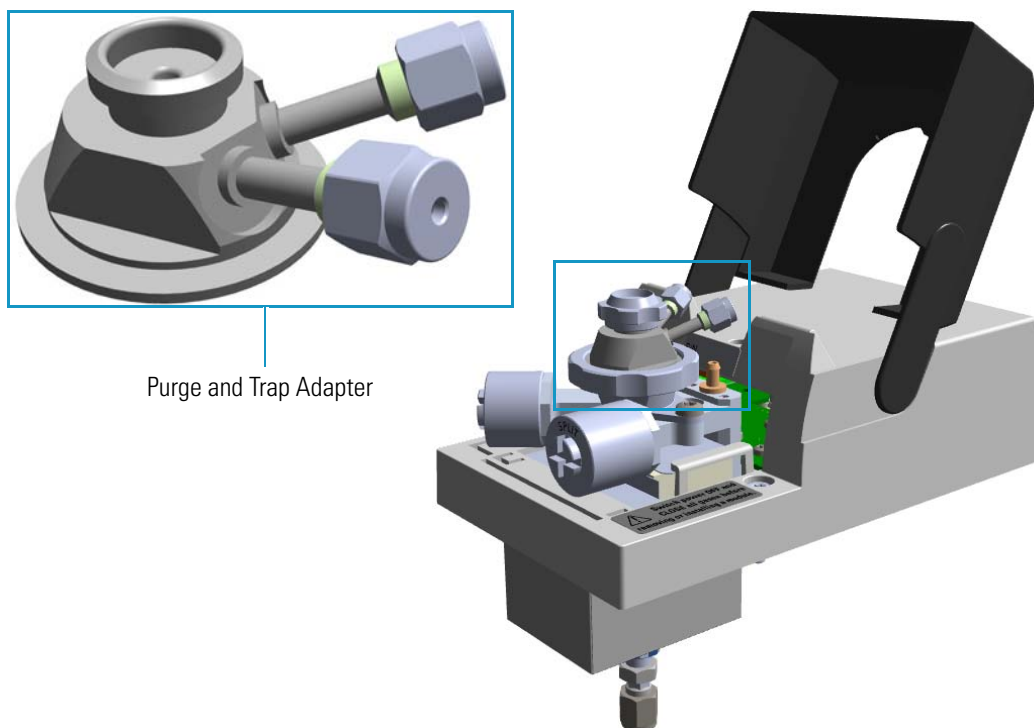

### ❖ To install the purge & trap adapter


1. Put the GC in standby condition.
2. Cool the oven, injector, and detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. Remove the column end from the injector.
  - a. Open the front door of the GC.


- b. Loosen the retaining nut from the injector fitting on the upper interior wall of the GC oven.
  - c. Remove the analytical column with its nut and ferrule from the bottom of the injector.
6. Put the autosampler away if present.
7. Remove the SSL/SSLBKF injector module from its seat.
  - a. Open the module flap cover.
  - b. Using a T20 Torxhead screwdriver, unscrew the three captive fixing screws.
  - c. Throw upward the module from its seat of the injector housing. Place the SSL/SSLBKF module on a clean surface.

**WARNING** Make sure the o-ring remains into its seat on the gas connection.  
**Do not install the module if the o-ring is missing.**

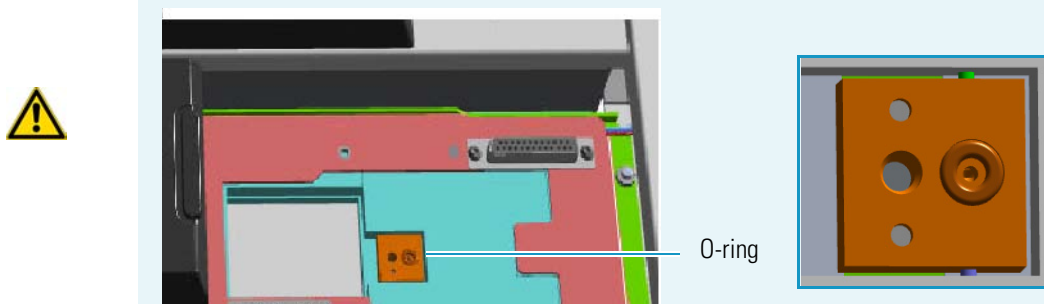

8. Replace the module flap cover with the one's provided. See [Figure 331](#).

**Figure 331.** Module Flap Cover for Purge & Trap

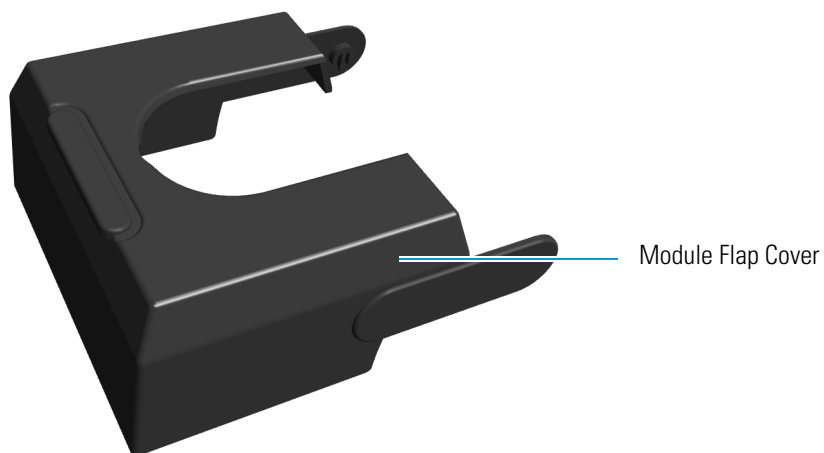

9. Reinstall the injector module into the main frame.
  - a. Open the module flap cover.

## 8 Installing Optional Kits

Installing the Purge & Trap Adapter Kit on the SSL/SSLBKF Injector


- b. Place the injector module in its seat. Be sure to insert the 25-pin male connector, on the bottom of the module, into the 25-pin female connector on the injector seat of the injector housing.
- c. Use a T20 Torxhead screwdriver to tighten the three captive fixing screws without overtighten.



**CAUTION** To maintain the correct alignment the screws must be tightened in turn. Tighten each screw only a small amount before moving to the next screw. Repeat until all are secure.

10. Remove the top parts of the injector. See [Figure 332](#).

**Figure 332.** Injector Top Parts Removal

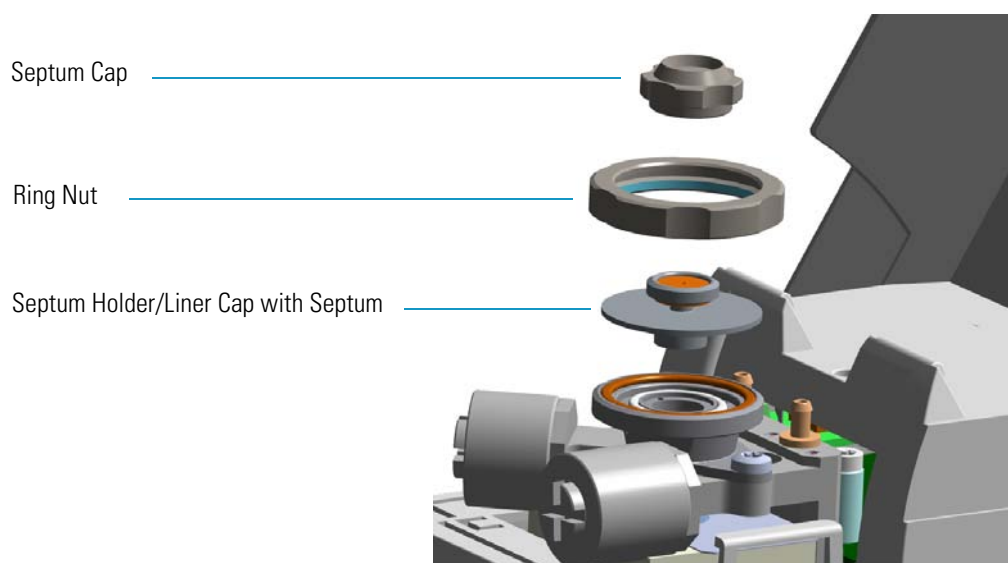

- a. Unscrew the septum cap of the injector.
- b. Unscrew the ring nut.
- c. Remove the septum holder/liner cap with septum from the injector body head.

11. Replace the liner.

**IMPORTANT** We suggest to replace the liner currently installed into the injector with the HS/SPME liner (PN 453A1335).


- a. Use tweezers to remove the current liner with the liner seal (o-ring) from the injector.

**CAUTION** Be careful not to break the glass liner when removing it. Glass splinters might fall into the lower part of the vaporization chamber.

- b. Holding the HS/SPME liner with tweezers place a new liner seal over the liner.
- c. Insert the liner into the injector and push it gently towards the bottom of the injector.


12. Install the purge & trap adapter on the top of the injector. See [Figure 333](#).

**Figure 333.** Purge & Trap Adapter (1)

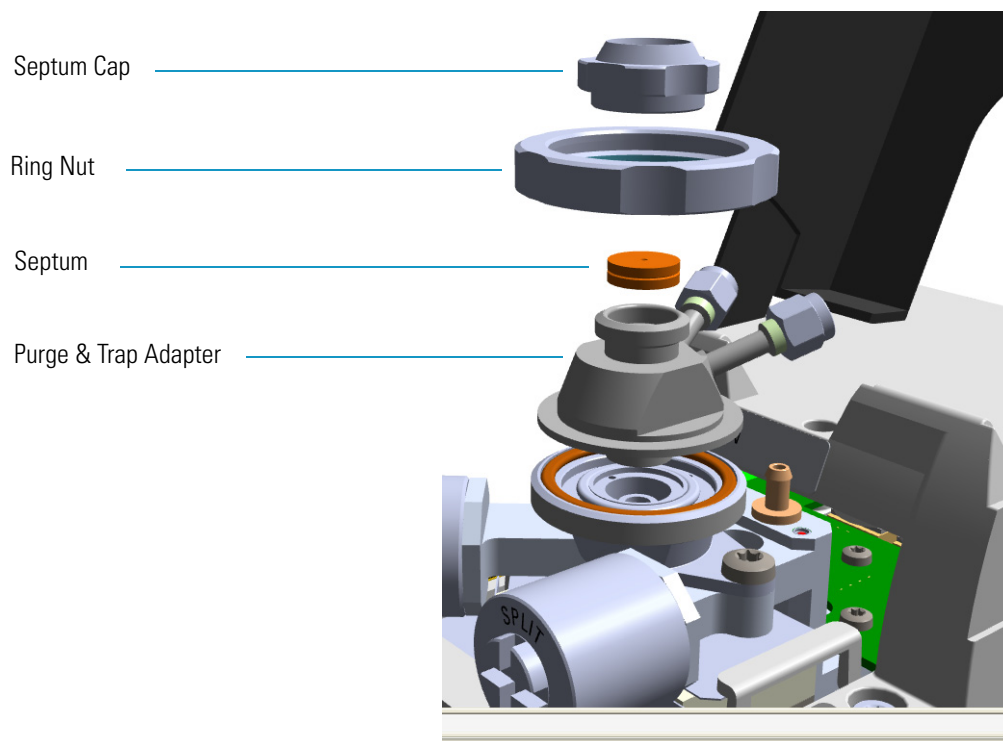

- a. Avoid touching the septum with your fingers. Insert a new septum into the septum holder cavity of the purge & trap adapter using tweezers. See [Figure 334](#).

**Figure 334.** Purge & Trap Adapter (2)

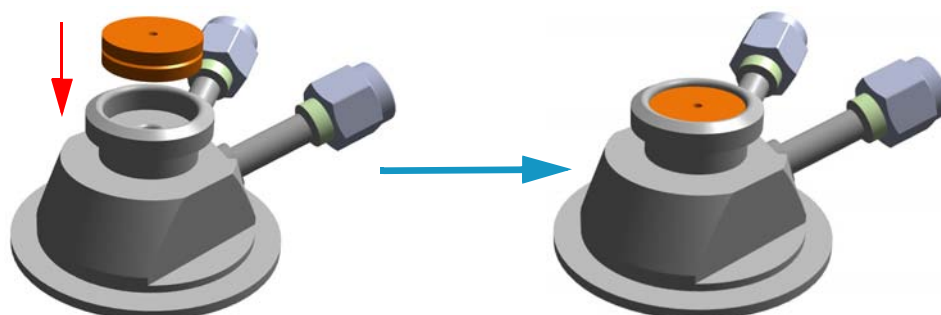


- b. Place purge & trap adapter with the septum on the body head of the injector. See [Figure 335](#).

## 8 Installing Optional Kits


Installing the Purge & Trap Adapter Kit on the SSL/SSLBKF Injector

**Figure 335.** Purge & Trap Adapter (3)

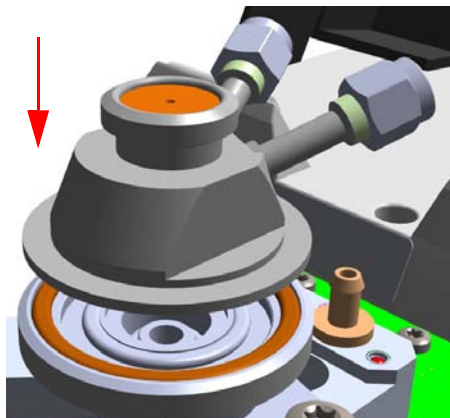


- c. Guide the ring nut on the purge & trap adapter with the septum, then fix it screwing the ring nut. See [Figure 336](#).

**Figure 336.** Purge & Trap Adapter (4)

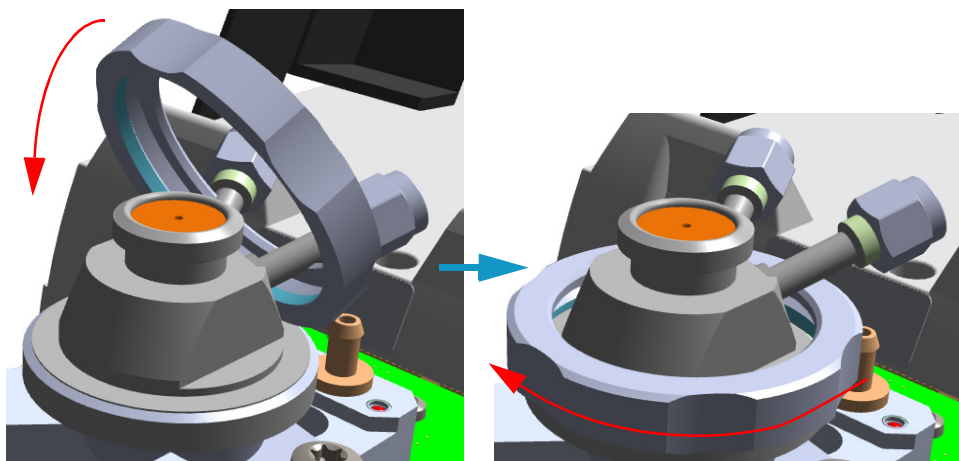

- d. Screw and tighten the septum cap to finger-tight. See [Figure 337](#).

**Figure 337.** Purge & Trap Adapter (5)

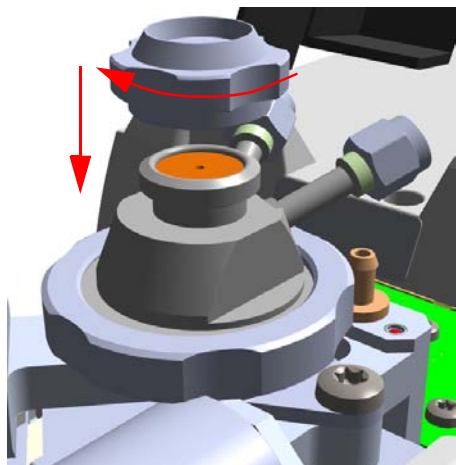




**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

13. Close the module flap cover.
14. Connect the purge & trap adapter to the Purge & Trap system. See [Figure 338](#).

**Figure 338.** Purge & Trap Adapter (6)

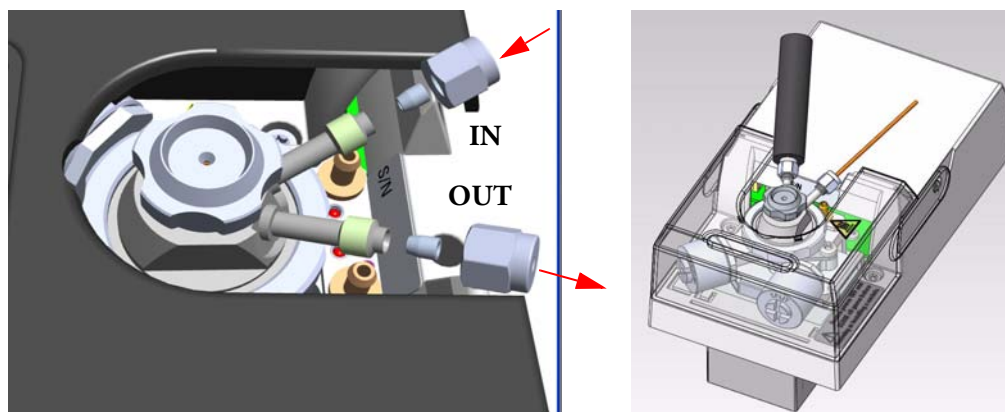

- a. Using the nut and ferrule provided, connect the heated transfer line from the Purge & Trap on the connector IN of the purge & trap adapter.
- b. Using the nut and ferrule provided, connect the gas line from the injector module to the Purge & Trap on the connector OUT of the purge & trap adapter.
15. Reconnect the column end to the injector and verify the connection point.
16. Open the gas supplies.
17. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip out the power switch (breaker) to the position I.
18. Pressurize the module with the carrier gas.
19. Check for leaks.
  - a. Use a handheld electronic leak detector (Thermo Scientific GLD Pro leak detector or equivalent) to check the two fittings for leaks.
  - b. If you detect a leak, tighten the connection and retest it.
  - c. Repeat this process until all connections are leak free.
20. Close the front door of the GC.
21. If present, update the autosampler for the new injection position.

## Installing the Packed Column Adapters

This section provides instructions for installing the adapters for connecting 1/8 or 1/16-in. OD metal packed columns to the bottom of the SSL injector, FID, NPD, ECD, TCD, and FPD. For the connection to a PDD detector see the **Caution** note below.



**WARNING** Do not try to connect packed columns if the lengths of the column ends above the end of the ferrule are not compatible with the adapters. In this case we suggest to connect a new packed column by using the nuts and ferrules provided with the kit.








**CAUTION Connecting a packed column to a PDD detector:** An adapter to connect the PDD packed column adapter can be assembled using a piece of tube, the 1/16-in. nut and gold plated ferrule, and a 1/8-in. to 1/16-in. reducing union in the standard outfit.

## Introduction

The adapters are provided with the kit PN 19050758. The kit includes the following parts:

**Table 14.** Packed Column Adapters and 1/2-in. Wrench
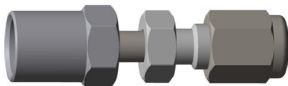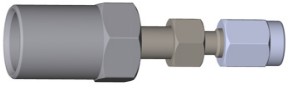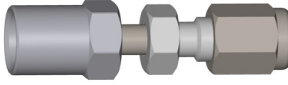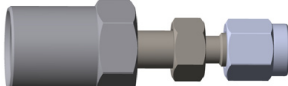
| Part  | Description   |
|---|---|
|  | Adapter for the connection of the 1/8-in. OD packed column to the SSL injector                        |
|  | Adapter for the connection of the 1/16-in. OD packed column to the SSL injector                       |
|  | Adapter for the connection of the 1/8-in. OD packed column to the FID, NPD, ECD, TCD, or FPD detector |
|  | Adapter for the connection of the 1/16-in. OD packed column to the FID, NPD, ECD, or FPD detector     |
|  | 1/2-in. wrench  |

## Getting Started

Install the column adapters following the instructions in the next operating procedures.

- “To install the packed column adapters” on [page 339](#)
- “To connect a new packed column to the injector and detector adapters” on [page 341](#)




❖ **To install the packed column adapters**

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Remove the analytical column.
  - a. Open the front door of the GC.
  - b. Loosen the retaining nut from the injector and detector fitting on the upper interior wall of the GC oven.
  - c. Remove the analytical column with its nut and ferrule from the bottom of the injector and the detector.
5. Remove the bottom parts of the SSL injector. See [Figure 339](#).

**Figure 339.** SSL Injector Bottom Parts Removal

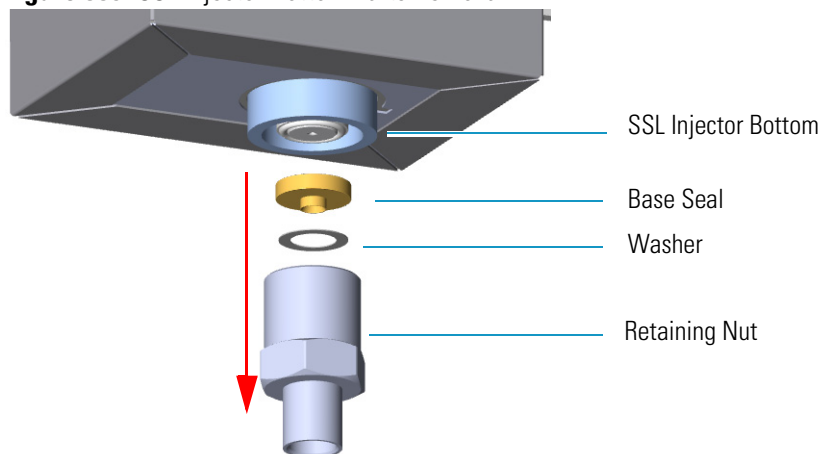


- a. Using the 1/2 wrench provided, unscrew the retaining nut with the base seal and the washer from the bottom of the injector.

**CAUTION** Make sure that the liner does not come out from the bottom of the injector.

Save the bottom parts of the SSL injector in a safe place because will be reused when you restore the original configuration for capillary columns.

6. Install the adapter on the bottom the injector interposing the washer. See [Figure 340](#).

**Figure 340.** Adapter Installation on SSL Injector Bottom

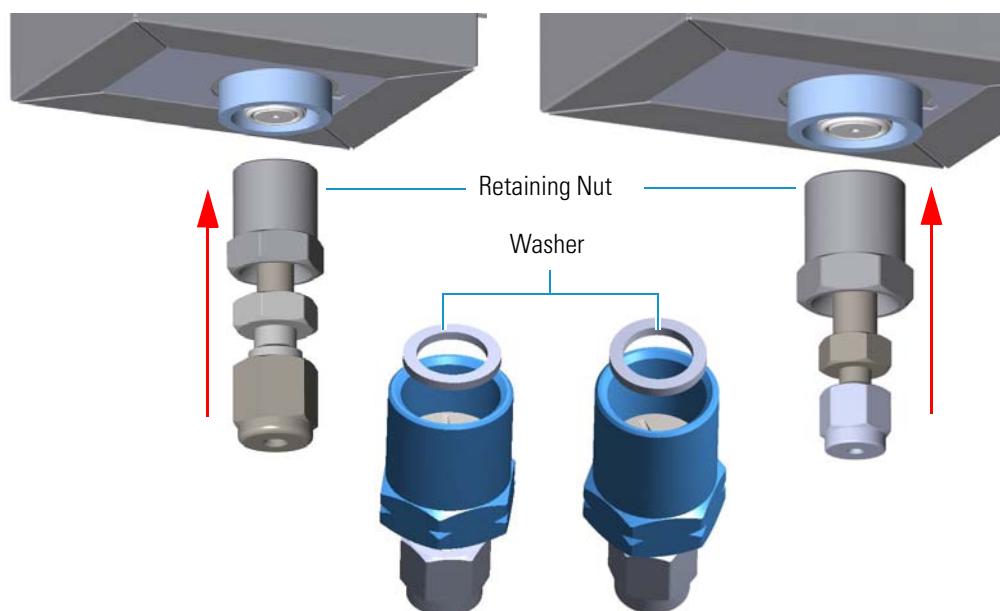


Adapter for 1/8-in. OD Packed Column

Adapter for 1/16-in. OD Packed Column

- a. Finger-tighten the retaining nut of the adapter until it start to grip the bottom of the injector.
  - b. Use the 1/2 wrench to firmly tighten the retaining nut of the adapter with the base seal from the bottom of the injector.
7. Remove the bottom parts of the FID, NPD, ECD, TCD, or FPD detector.  
See [Figure 341](#).

**Figure 341.** Detector Bottom Parts Removal

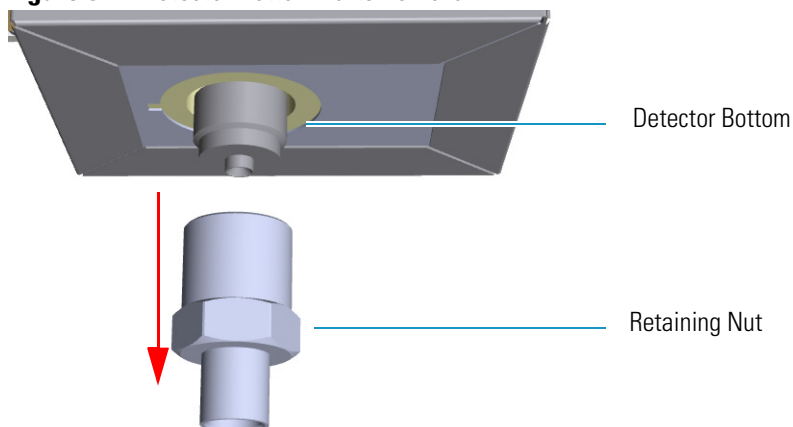


**Note** Save the bottom parts of the detector in a safe place because will be reused when you restore the original configuration for capillary columns.

- a. Using the 1/2 wrench provided, unscrew the retaining nut from the bottom of the detector.

8. Install the adapter on the bottom the detector interposing the washer. See [Figure 342](#).

**Figure 342.** Adapter Installation on the Bottom of the Detector

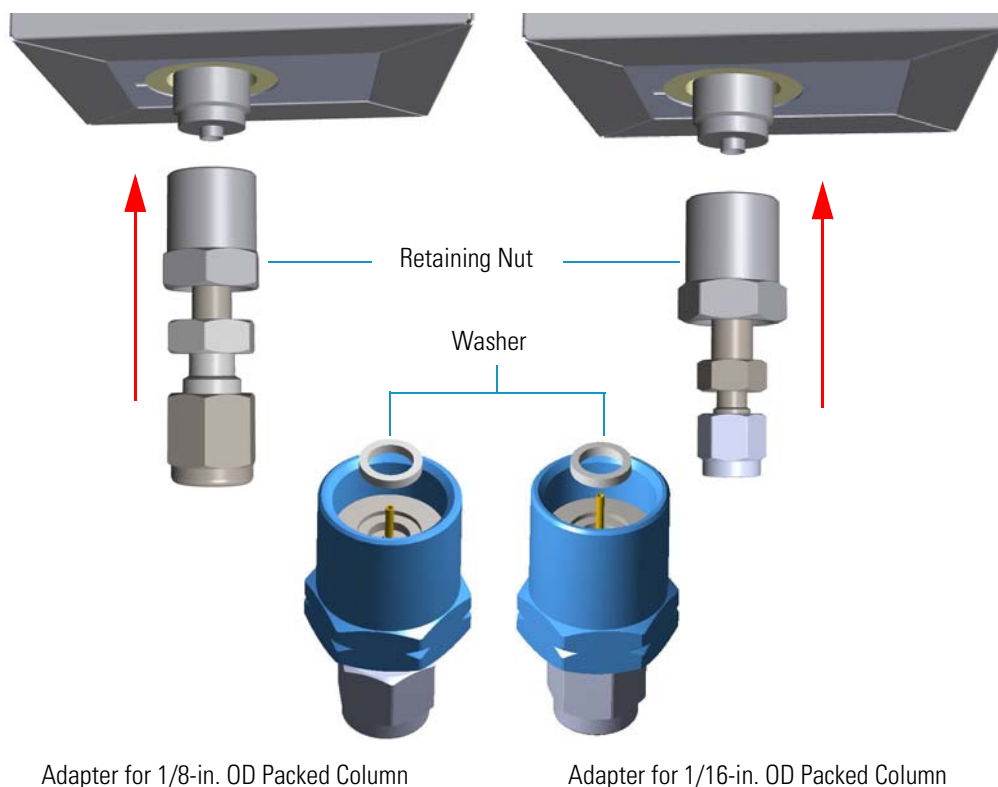


- a. Finger-tighten the retaining nut of the adapter until it start to grip the bottom of the detector.
- b. Use the 1/2 wrench to firmly tighten the retaining nut of the adapter with the base seal from the bottom of the injector.

❖ **To connect a new packed column to the injector and detector adapters**


Before you begin, verify that the proper adapters are installed on the injector and detector side.

1. Preparing the metal packed column.
  - a. Using a 1/4-5/16-in. wrench, remove the column nut and ferrules from the base of the injector adapter.
  - b. Using a 1/4-5/16-in. wrench, remove the column nut and ferrules from the base of the detector adapter.
  - c. Slide the column nut and the ferrules onto the packed column injector and detector ends in the order and direction as shown in the [Figure 343](#) and [Figure 344](#).

**Figure 343. Column Nut and Ferrules (1)**

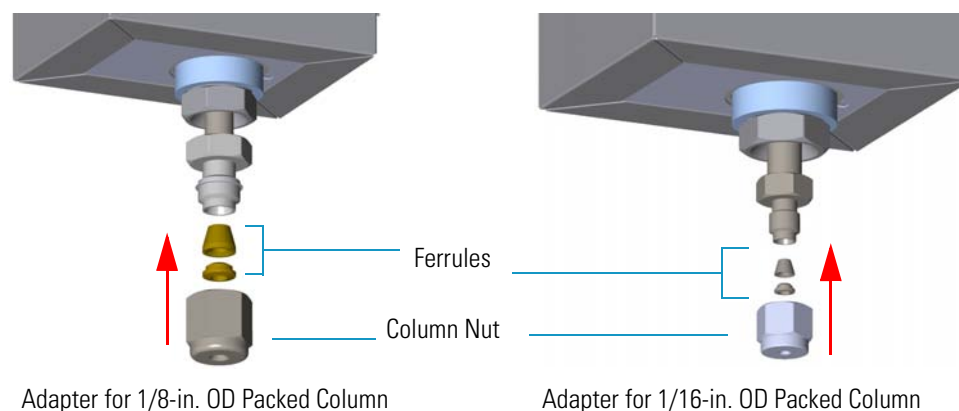

**Figure 344. Column Nut and Ferrules (2)**

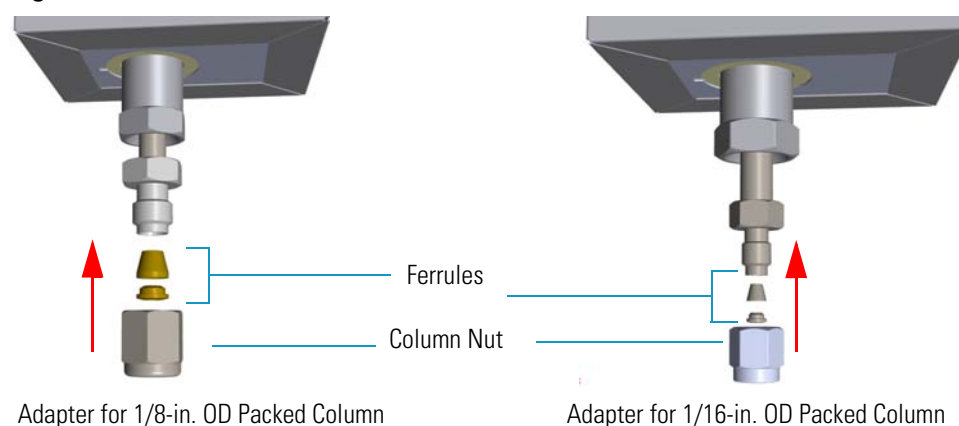

2. Connect the packed column to the injector.
  - a. Insert the inlet end of the column as far as possible into bottom of the adapter.
  - b. Slide the ferrule up to adapter base then finger-tighten the column retaining nut until it starts to grip the column.
  - c. Finger-tighten the retaining nut until it starts to grip the column plus a quarter turn.
3. Close the front door of the GC.
4. Turn on the gas supplies and the injector temperature.
5. Setup the GC.
  - a. Set the oven and injector temperature to 50 °C.
  - b. Allow the column to purge for few minutes.
6. Condition the column.

The column must be conditioned before inserting it into the detector.

Column conditioning consists of passing a carrier gas flow through the column heated at a programmed temperatures as described in the *column manufacturer's instructions*.

In case the column does not have any column conditioning instructions, perform the column conditioning by setting a final temperature up to 10 °C - 20 °C below its recommended maximum temperature.



**CAUTION** When performing column conditioning, the column should be connected only to the injector leaving the column outlet disconnected to avoid the possibility of contamination of the detector.

Do not use hydrogen as the carrier for conditioning! It could vent into the oven and present an explosion hazard.

- a. Run the temperature program that is recommended by the manufacturer.



**INSTRUMENT DAMAGE:** Never exceed the column manufacturer's maximum operating temperature.

7. Connect the column to the detector inside the GC.

- a. Lower the oven temperature to 30 °C and allow it to cool.



**WARNING-BURN HAZARD:** The injector, detector, and oven, may be hot. Allow them to cool to room temperature before touching them.

8. Connect the packed column to the detector.

- a. Insert the inlet end of the column as far as possible into bottom of the adapter.
- b. Slide the ferrule up to adapter base then finger-tighten the column retaining nut until it starts to grip the column.
- c. Use the wrench to tighten the retaining nut. Use no more pressure than is necessary to obtain a good seal (1/4 to 1/2 turn).


9. End of the column installation.

- a. Close the front door of the GC.

# Installing the HS Adapter Kit on the SSL/SSLBKF Injector

This section provides the instructions for installing the HS Adapter kit PN 19050732, on the SSL/SSLBKF injector for the use with a TriPlus 300 Headspace sampling system. See [Figure 345](#).

**Figure 345.** HS Adapter

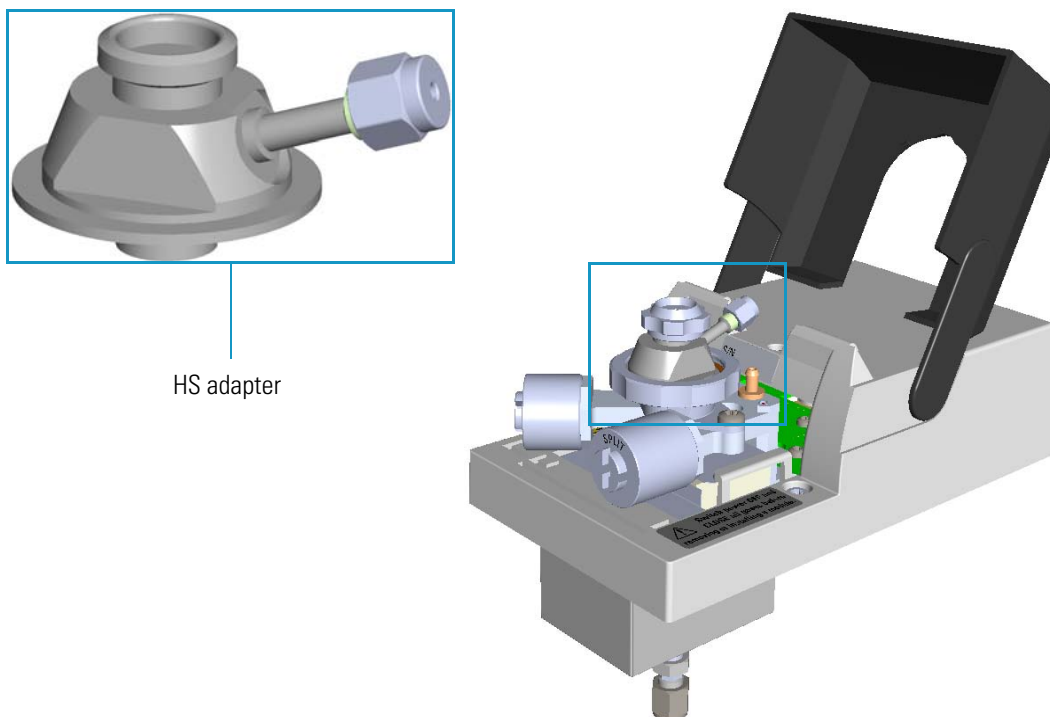

### ❖ To install the HS adapter


1. Put the GC in standby condition.
2. Cool the oven, injector, and detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. Remove the column end from the injector and the detector.
  - a. Open the front door of the GC.


- b. Loosen the retaining nut from the injector fitting on the upper interior wall of the GC oven.
  - c. Remove the analytical column with its nut and ferrule from the bottom of the injector.
6. Remove the SSL/SSLBKF injector module from its seat.
  - a. Open the module flap cover.
  - b. Using a T20 Torxhead screwdriver, unscrew the three captive fixing screws.
  - c. Lift the module from its seat in the injector housing. Place the SSL/SSLBKF module on a clean surface.

**WARNING** Make sure the o-ring remains into its seat on the gas connection.  
**Do not install the module if the o-ring is missing.**

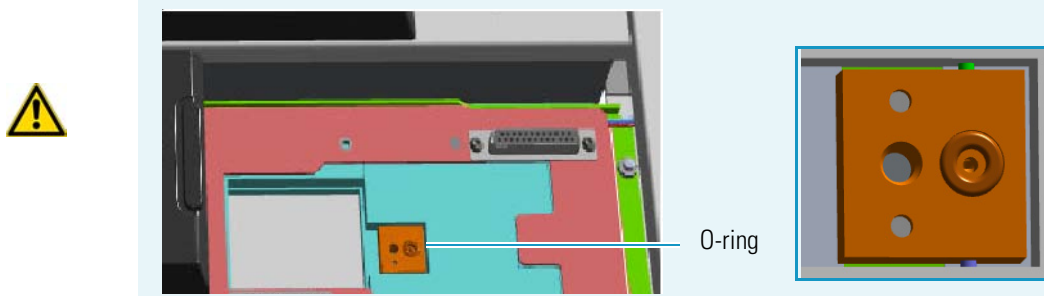

7. Replace the module flap cover with the one's provided. See [Figure 346](#).

**Figure 346.** Module Flap Cover for the HS Adapter

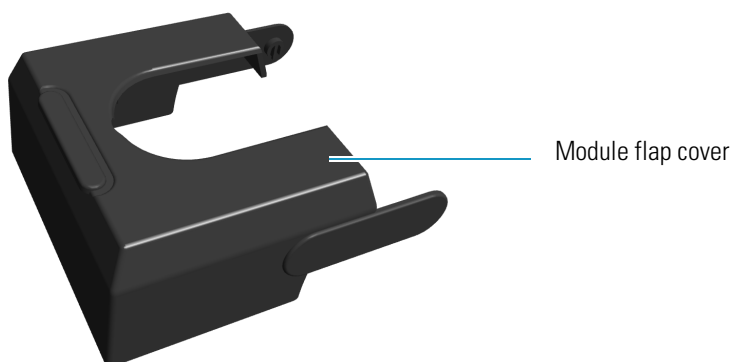

8. Reinstall the injector module into the main frame.
  - a. Open the module flap cover.
  - b. Place the injector module in its seat. Be sure to insert the 25-pin male connector, on the bottom of the module, into the 25-pin female connector on the injector seat of the injector housing.
  - c. Use a T20 Torxhead screwdriver to tighten the three captive fixing screws without overtightening.

## 8 Installing Optional Kits


Installing the HS Adapter Kit on the SSL/SSLBKF Injector



**CAUTION** To maintain the correct alignment the screws must be tightened in turn. Tighten each screw only a small amount before moving to the next screw. Repeat until all are secure.

9. Remove the top parts of the injector. See [Figure 347](#).

**Figure 347.** Injector Top Parts Removal

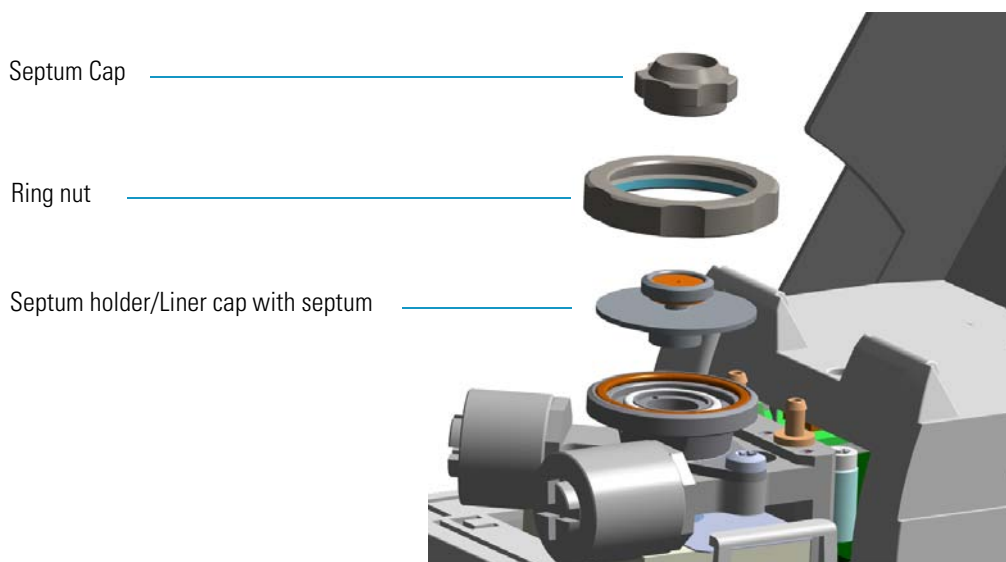

- a. Unscrew the septum cap of the injector.
- b. Unscrew the ring nut.
- c. Remove the septum holder/liner cap with septum from the injector body head.

10. Replace the liner.

**IMPORTANT** We suggest replacing the liner currently installed into the injector with the HS/SPME straight glass empty liner PN 453A1335.


- a. Use tweezers to remove the current liner with the liner seal (o-ring) from the injector.

**CAUTION** Be careful not to break the glass liner when removing it. Glass splinters might fall into the lower part of the vaporization chamber.


- b. Holding the HS/SPME liner with tweezers, place a new liner seal over the liner.
- c. Insert the liner into the injector and push it gently towards the bottom of the injector.

11. Install the HS adapter on the top of the injector. See [Figure 348](#).

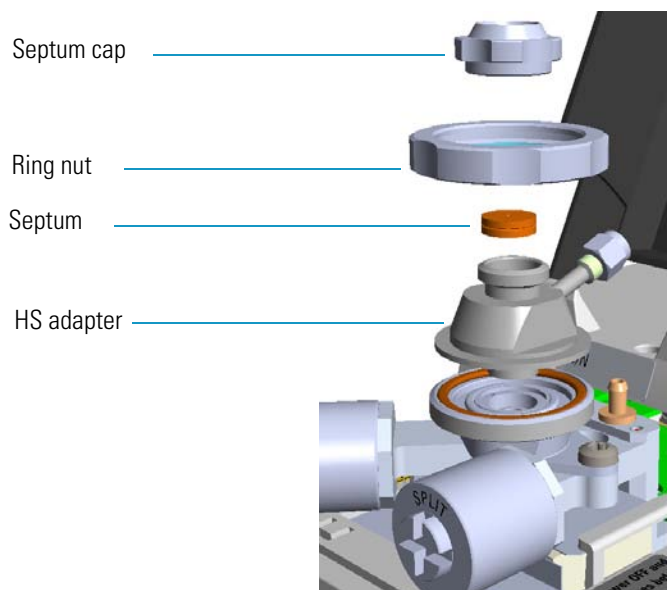


**Figure 348.** HS Adapter (1)

- a. Avoid touching the septum with your fingers. Insert a new septum into the septum holder cavity of the HS adapter using tweezers. See [Figure 334](#).

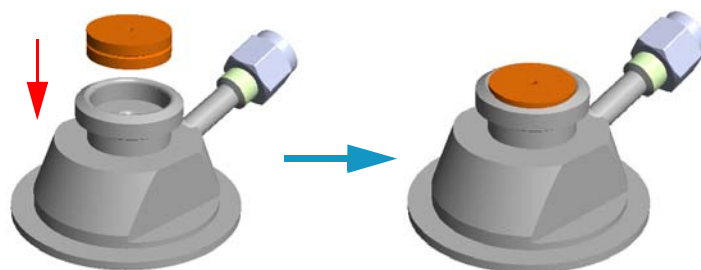**Figure 349.** HS Adapter (2)

- b. Place the HS adapter with the septum on the body head of the injector. See [Figure 350](#).

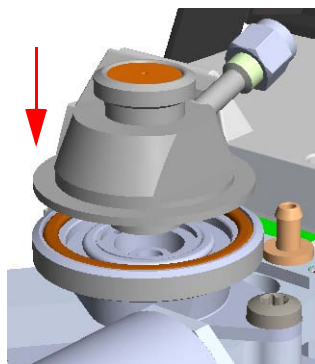
**Figure 350.** HS Adapter (3)

- c. Guide the ring nut on the HS adapter with the septum, then fix it by screwing in the ring nut. See [Figure 351](#).

## 8 Installing Optional Kits


Installing the HS Adapter Kit on the SSL/SSLBKF Injector

**Figure 351.** HS Adapter (4)

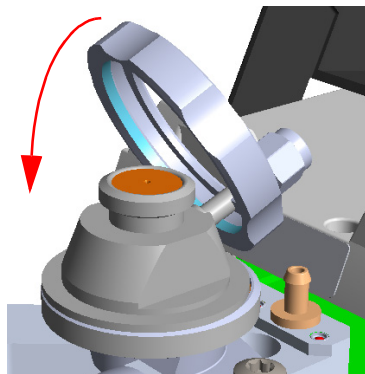


12. Close the module flap cover.
13. Connect the carrier gas line coming from the TriPlus 300 Headspace sampling system to the HS adapter by using the proper nut and ferrule. See [Figure 352](#).

**Figure 352.** HS Adapter (5)

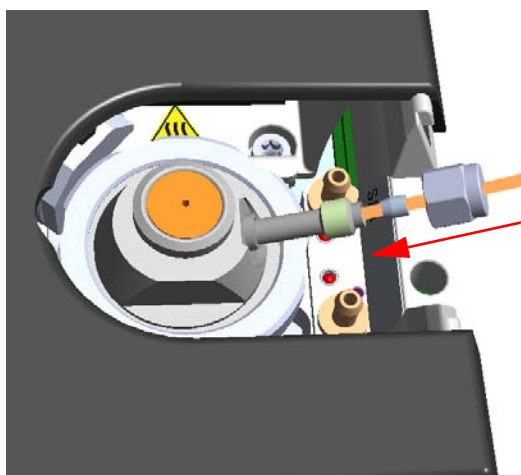

14. Screw and finger-tighten the septum cap provided. See [Figure 353](#)

**Figure 353.** HS Adapter (6)

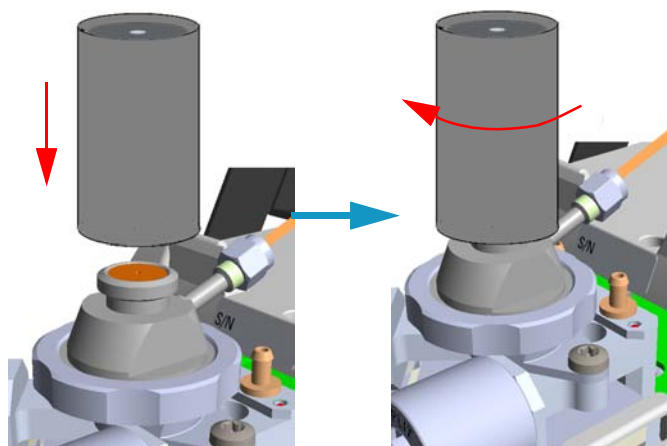




**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

15. Insert the transfer line needle into the septum cap. See [Figure 354](#).


**Figure 354.** Transfer Line Connection

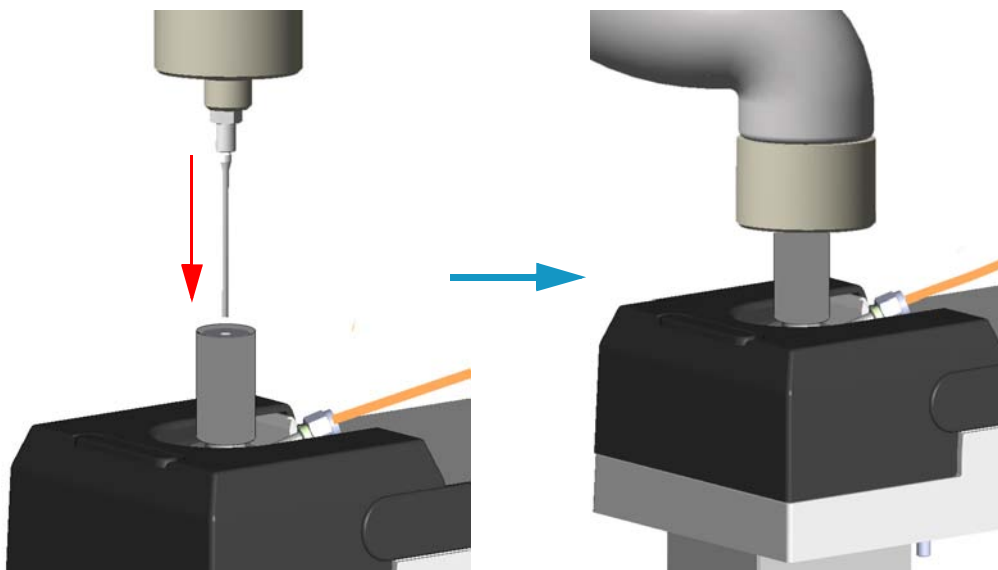

16. Cut the piece of column end previously connected to the injector.
- Use a scoring wafer to score and break the column in order to remove the current ferrule and the nut.
17. Connect the column to the injector.
- Wipe about 100 mm (4 in.) of the column with a tissue soaked in methanol.
  - Insert the column through the injector retaining nut and the proper ferrule (open end up). Wipe the column again with a tissue soaked in methanol.
  - Use a scoring wafer to score and break the column about 1 cm (0.4 in.) from the end. Use a magnifying glass to check for an even, flat cut. Repeat if necessary.
- Tip** Slide a notched septum on the column before the injector retaining nut to make it easier to measure the proper distance between the nut and the end of the column.
- Position the column so that the end of the column extends the proper distance of **30 mm** above the end of the ferrule.
  - Insert the notched septum on the column to hold the retaining nut at this position. Thread the retaining nut into the injector but do not tighten.

## Installing the HS Adapter Kit on the SSL/SSLBKF Injector

- Figure 355.** Using a Flowmeter for Leak Check

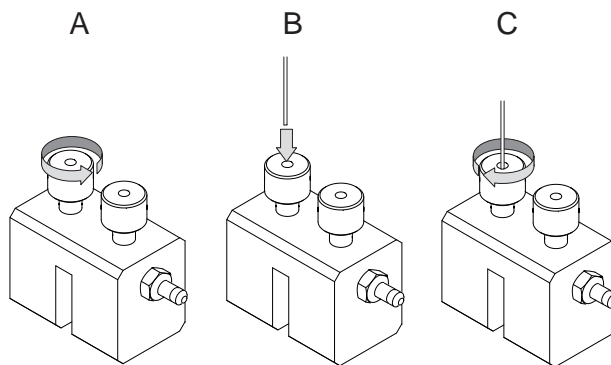

- b. If your GC is equipped with the touch screen as user interface, select the **Leak Check** icon in the **Maintenance** menu, otherwise perform the Leak Check through the Chromatography Data System by selecting the proper function.
- c. Start the leak check to begin operation. The split and purge valves of the selected channel are automatically closed and the channel is pressurized with carrier gas to the leak check setpoint.


- d. The system monitors the pressure for one minute. If the pressure does not drop more than the maximum allowed sensitivity value, then the leak check will pass.

If the leak check did not pass, you should use the leak detector to find and fix the leaks.

**Tip** Leaks can be caused by not tightening the fitting on the column-flowmeter connector. We recommend that you check that fitting before looking elsewhere.

- e. Repeat the leak check until no leaks are indicated.
22. Calibrate the carrier gas flow (column evaluation).
    - a. Carefully push the capillary column end into the flow meter section of the column-flowmeter connector. See [Figure 356](#).

**Figure 356.** Column Flowmeter Connector

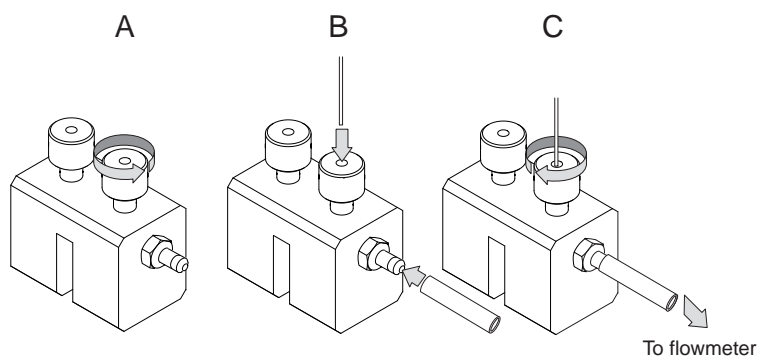

- b. Connect the flowmeter to the dedicated fitting on the column-flowmeter connector.
- c. If your GC is equipped with the touch screen as user interface, select **Back** or **Front Column** icon in the **Configuration** menu, otherwise perform the Column Evaluation through the Chromatography Data System by selecting the proper function.
- d. Select the column and input the physical characteristics of the column.
- e. If a pre-/post-column is present, set the length and nominal internal diameter of the pre-/post-column in the same valid ranges for the column. The following two lines are added to the menu.
- f. According to the physical characteristics of the column, the system calculates and displays the relevant Column K-factor.

**Note** For the most reproducible results, you should conduct a more detailed column evaluation. However, the following steps, while recommended, are not required.

- g. Start column evaluation. At the end of the routine, a message will indicate that evaluation was successful.

## 8 Installing Optional Kits

Installing the HS Adapter Kit on the SSL/SSLBKF Injector

- h. Expect a K-factor of approximately 0.7 – 0.9 for a 15 m, 0.25 mm ID column (1.3 – 2.0 for a 30 m, 0.25 mm ID column). If the column does not report a K-factor within this range or within 0.1 units of the previous stored value, check for a leak or broken column using the leak detector. The K-factor is a measured resistance for the column. A K-factor that is too low may indicate a leak in the system, while a K-factor that is too high may indicate a blockage.

### 23. Disconnect the column-flowmeter.

- a. Disconnect the column from the column-flowmeter connector.
- b. Remove the clear plastic component, including its fittings, from the oven and set it aside.
- c. Close the GC door.

### 24. Condition the column.

The column must be conditioned before inserting it into the detector.

Column conditioning consists of passing a carrier gas flow through the column heated at a programmed temperatures as described in the *column manufacturer's instructions*.

In case the column does not have any column conditioning instructions, perform the column conditioning by setting a final temperature up to 10 °C - 20 °C below its recommended maximum temperature.



**CAUTION** When performing column conditioning, the column should be connected only to the injector leaving the column outlet disconnected to avoid the possibility of contamination of the detector.

Do not use hydrogen as the carrier for conditioning! It could vent into the oven and present an explosion hazard.

- a. Run the temperature program that is recommended by the manufacturer.



**INSTRUMENT DAMAGE:** Never exceed the column manufacturer's maximum operating temperature.

### 25. Connect the column to the detector inside the GC.

- a. Lower the oven temperature to 30 °C and allow it to cool.



**WARNING-BURN HAZARD:** The injector, detector, oven, and transfer line may be hot. Allow them to cool to room temperature before touching them.

- b. Unwind the column enough to easily connect its ends to the injector and the detector.

**Note** Wear clean, lint- and powder-free gloves when you handle the column and injector ferrule.

- c. Wipe about 100 mm (4 in.) of the column with a tissue soaked in methanol.
- d. Use a scoring wafer to score and break the column outlet about 2.5cm (1 in.) from the end. Use a magnifying glass to check for an even, flat cut. Repeat if necessary.
- e. Insert the column through the proper detector retaining nut and ferrule (open end up). Wipe the column again with a tissue soaked in methanol.

**Tip** Slide a notched septum on the column before the detector retaining nut to make it easier to measure the proper distance between the bottom nut and end of the column.

- f. Position the column so that the end of the column extends the proper distance above the end of the ferrule as reported in [Table 15](#).  
For **PDD** see the instruction described at the [step g](#) on [page 353](#).

**Table 15.** Column Insertion Depth For FID, NPD, TCD, ECD, FPD, and PDD Detectors

| FID  | NPD | TCD | ECD   | FPD    | PDD    |
|--|-----|-----|-------|--------|--------|
| Insert the column as far as goes and withdrawn about 2mm |     |     | 23 mm | 125 mm | 136 mm |

- i. For **FID**, **NPD**, **TCD**, insert the column into the detector, paying attention to not force it further. Finger-tighten the retaining nut, then withdraw the column 2-3 mm. Tighten the retaining nut an additional a quarter turn.
- ii. For **ECD** and **FPD**, insert the notched septum on the column to hold the retaining nut in this position. Thread the retaining nut into the detector but do not tighten.  
Finger-tighten the retaining nut until it starts to grip the column plus a quarter turn. Remove the notched septum from the column.
- g. For **PDD** the column must penetrate **136 mm** inside the capillary column adapter.

**Note** When inserting the capillary column into the PDD detector it might rarely happen to feel a slight resistance. In this case, for proper column installation, pull the column out slightly and adjust the angle before inserting it further.



**IMPORTANT** To install a packed column, the pre-installed capillary column adapter must be replaced with the **packed columns adapter** that enters into the PDD cell for the correct length.

26. End of the column installation.
  - a. Close the front door of the GC.

## Installing the Large Volume Splitless Kit

This section provides instruction for installing the Large Volume Splitless kit (PN 19050725) on your TRACE 1300/TRACE 1310.

The Large Volume Splitless injector is a setup of the standard splitless injector, where the introduction of large amount of liquid samples can be performed manually, or with the TriPlus RSH, TriPlus 100 Liquid Sampler, or AI/AS 1310 autosampler.


Large Volume-Splitless kit includes:

- Two dummy filters
- Two o-rings for dummy filters
- An deactivated connector (press-fit)
- An uncoated precolumn (5 m x 0.32 mm ID)
- A dedicated splitless liner (set of 5)
- LV-SL Assistant software


### ❖ To install the large volume splitless kit

1. Put the GC in standby condition.
2. Cool the oven and injector to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.
4. Put the autosampler away if present.
5. Open the module flap cover
6. Replace the filters with the dummy filters.
  - a. Remove both the filters from their seats by turning them counter-clockwise.



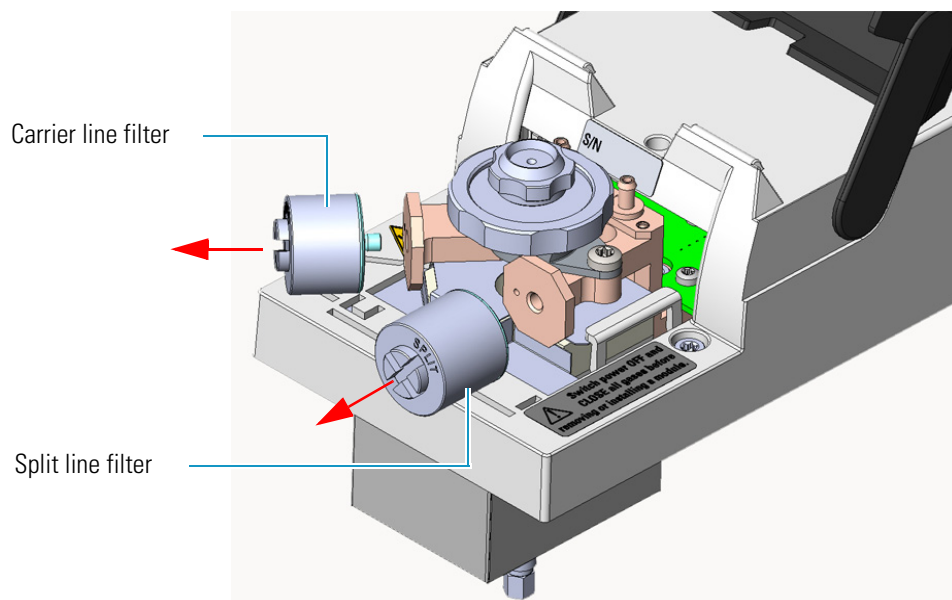


- b. Install the dummy filters in their seats interposing the o-ring, then turn them clockwise.

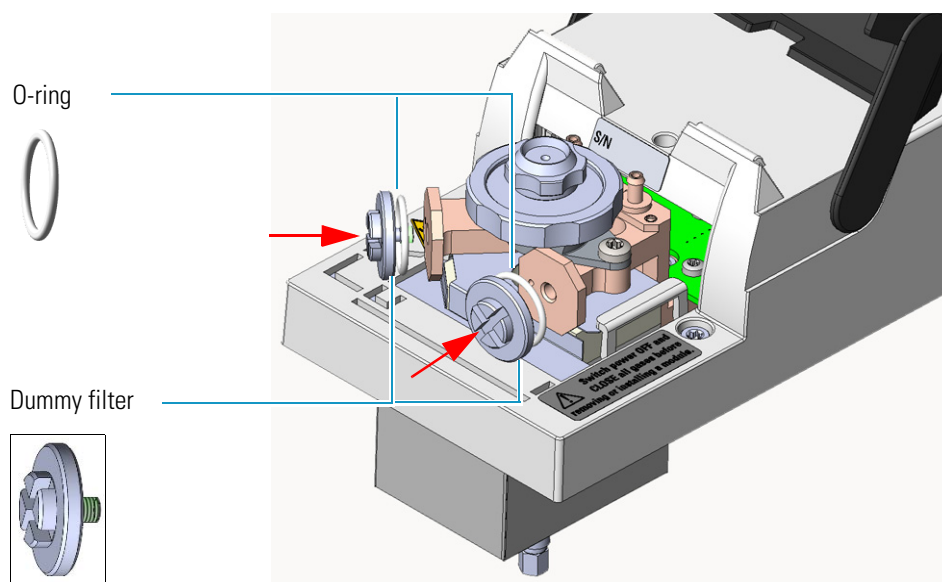

7. Replace the current liner installed into the Split/Splitless injector with the dedicated splitless liner.
8. Close the module flap cover.
9. Install the uncoated precolumn.
  - a. Open the oven door.
  - b. Disconnect the analytical column from the bottom of the injector.

- c. Connect an end of the precolumn to the bottom of the injector by using the proper nut and ferrule. Position the precolumn so that its end extends a distance of 5 mm above the end of the ferrule.


10. Couple the precolumn to the analytical capillary column.

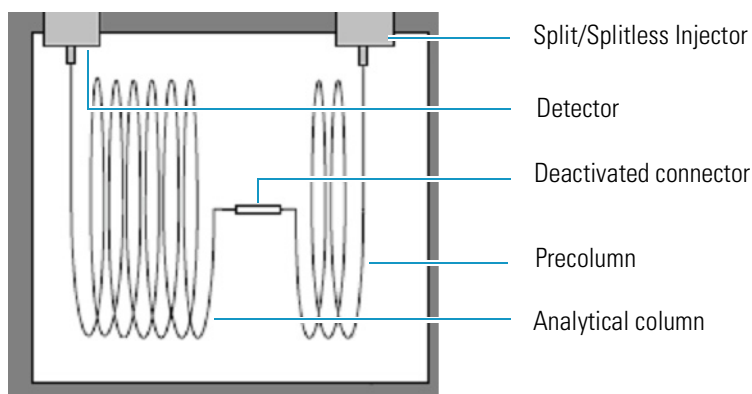

- a. Properly cut the fused silica column ends pay attention to achieve a clean square cut by using a ceramic scoring wafer or sapphire scribe.

**CAUTION** A poorly cut will produce an insufficient seal.

- b. insert the precolumn and analytical column ends into the relevant ports of the deactivated connector.

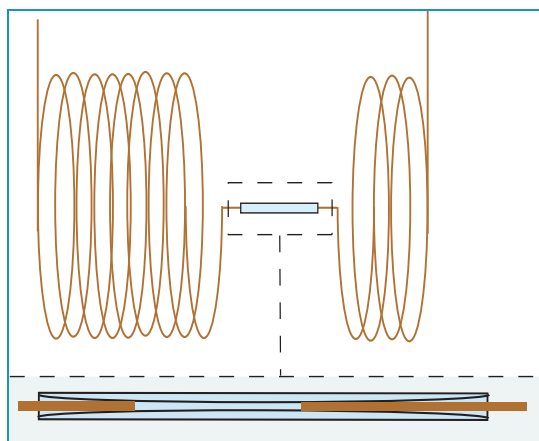


**Note** To create a good seal between all the parts, will be necessary to increase the oven temperature up to 200 °C.

11. Close the oven door.
12. If present, move the autosampler towards the module to restore the original alignment.
13. Turn the carrier gas on.
14. Set the injector, detector, and GC working conditions.

## Installing the Manual On/Off Valve for Single Gas Line

This section provides instructions for installing the manual on/off toggle valve for single gas line kit (PN 19050756) on a GC gas inlet of your TRACE 1300/TRACE 1310. See [Figure 357](#).

**Figure 357.** Manual Toggle Valve Installation

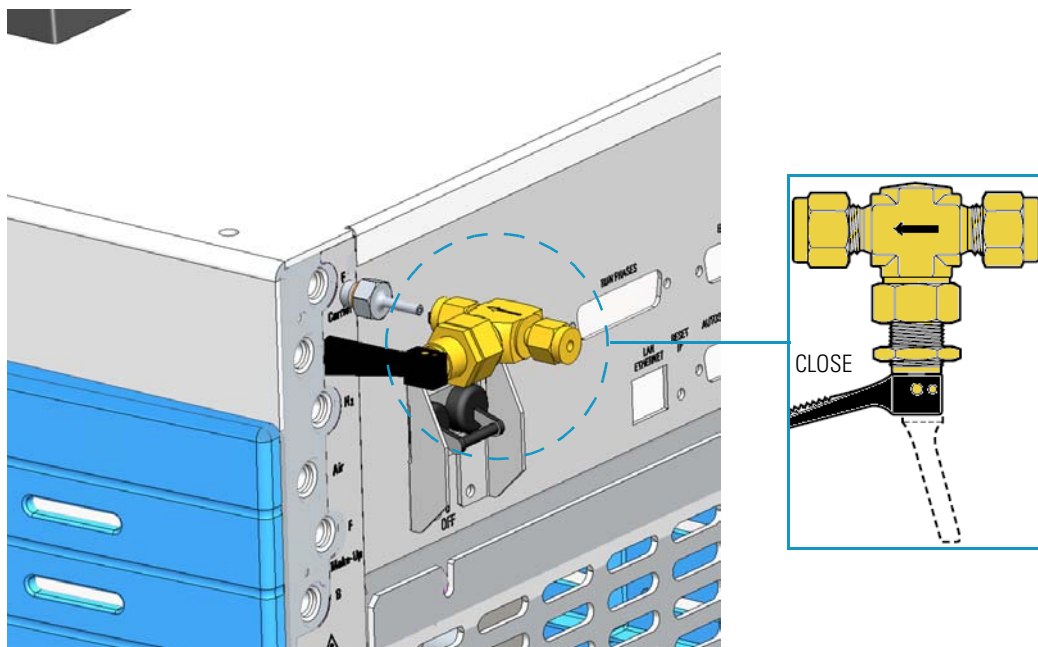

The manual toggle valve is used for manually opening and closing the gas flow coming from the supply line.

It is particularly useful when for any reason you need to manually interrupt the flow of a gas into the GC, for example: before the replacement of a module, to save the consumption of a gas when its flow is unnecessary, and so on.

The kit can be installed only on a single gas line, then if you need managing more gas lines, up to six kits can be installed accordingly.

### ❖ To install the manual on/off valve kit


**Note** This procedure describes the installation of a manual on/off valve on a single gas line. In the example the installation of the valve on the Front Carrier gas line is considered.

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.
3. Turn the carrier gas off, and wait for the carrier pressure to go to zero.
4. On the back of the GC disconnect the gas supply tube and the gas inlet fitting with its o-ring from the GC gas inlet where the manual on/off valve must be installed.


## 8 Installing Optional Kits

### Installing the Manual On/Off Valve for Single Gas Line


Use a 7/16-in. wrench for unscrewing the fittings.

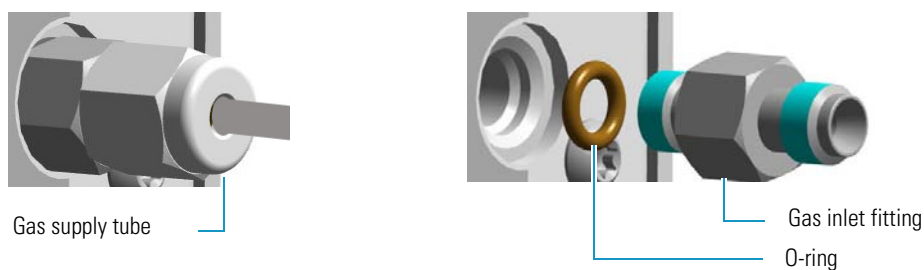


5. Screw the gas inlet fitting 1/8-in. provided into the gas inlet port interposing the o-ring. Use a 7/16-in. wrench for tightening the fitting.

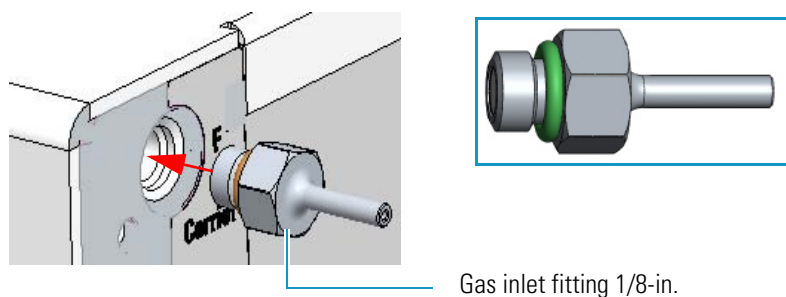


6. Using the nut and ferrules provided, connect the valve outlet to the gas inlet fitting 1/8-in. Make sure that the arrow marked on the valve body indicating the direction of the flow is turned forward the GC gas inlet. Use a 7/16-in. wrench for tightening the fittings.

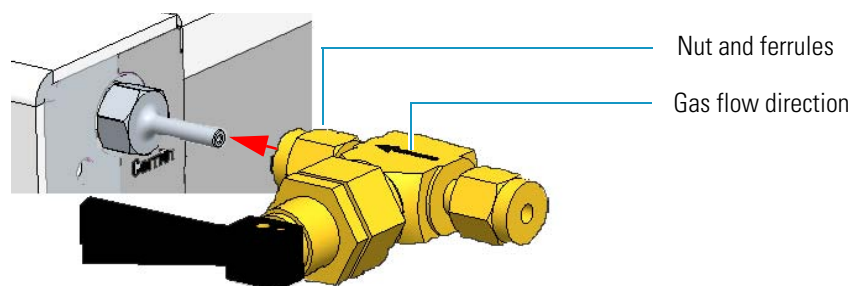

7. Remove the nut and the ferrules provided from the valve inlet.

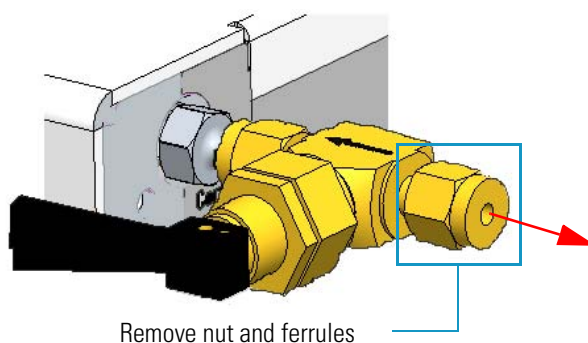

8. Connect the gas supply tube previously disconnected (see the step 4). Use a 7/16-in. wrench for tightening the fittings.

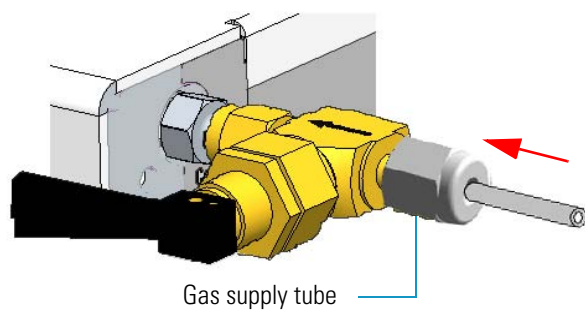

9. Turn the carrier gas on.
10. Open or close the flow of the supply gas turning the black handle of the valve up or down respectively.
11. Set the injectors, detectors, and GC working conditions.

## Connecting a SSL/PTV Backflush System

This section provides instructions for connecting your SSLBKF/PTVBKF injector modules with precolumn and analytical capillary column into the oven. Refer to the chapters SSLBKF Injector Module and PTVBKF Injector Module in the *TRACE 1300/TRACE 1310 User Guide*. This system operates up to 300 °C.


The operation consists of the following operating procedures:

- “To assemble the Tee connector on the mounting bracket” on page 360
- “To install the mounting bracket” on page 367
- “To connect backflush line, precolumn, and capillary column” on page 369
- “To install of a second backflush system” on page 373




**CAUTION** Before starting make sure the SSLBKF/PTVBKF injector module is correctly installed into its seat.

### ❖ To assemble the Tee connector on the mounting bracket

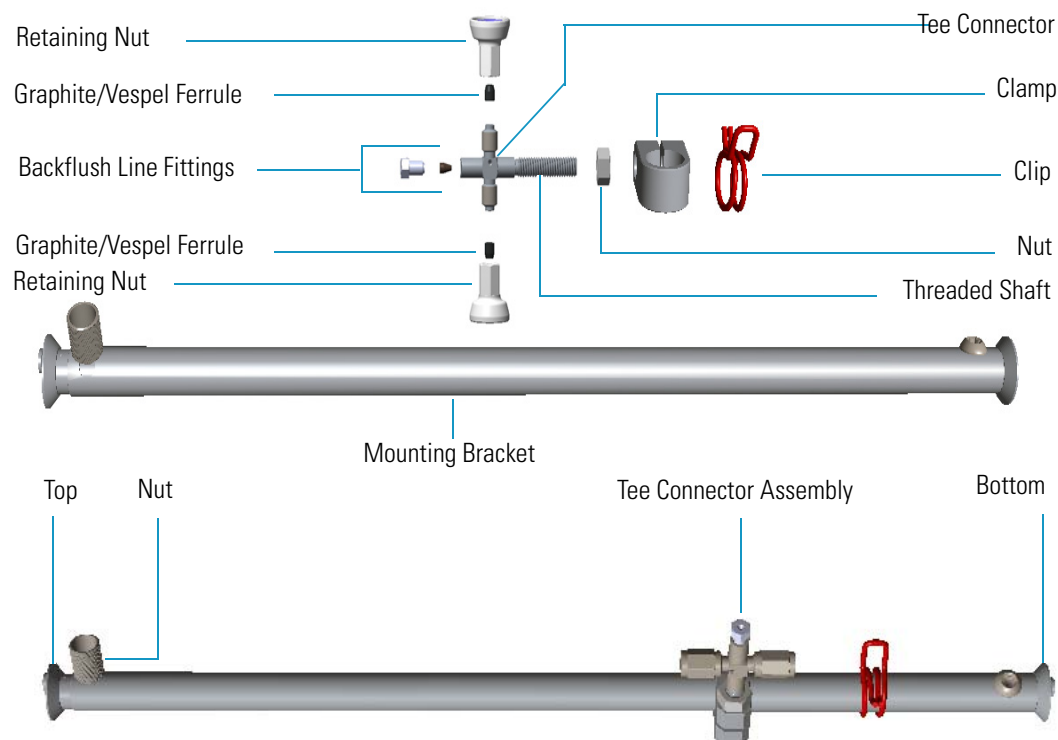

1. Remove the nut and the top end of the mounting bracket. This nut secures the two sections of the mounting bracket. See [Figure 358](#).


**Figure 358.** Assembling Tee Connector and Mounting Bracket (1)

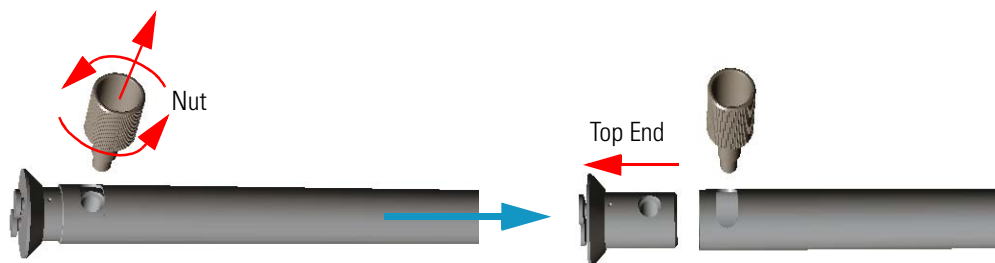

**Note** Wait reinstalling the top end of the mounting bracket and the nut. These parts will be reinstalled during the installation of the mounting bracket into the GC oven. See [“To install the mounting bracket”](#) on page 367.


- Slide the clip and the clamp on the body of the mounting bracket. See [Figure 359](#).

**Figure 359.** Assembling Tee Connector and Mounting Bracket (2)

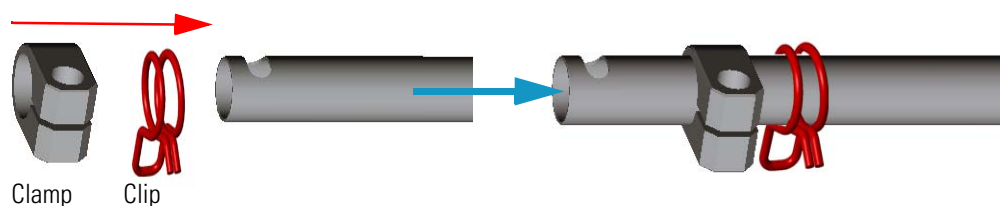


- Screw the part of threaded shaft of the Tee connector with the smaller diameter into the Tee connector, then screw the nut on the threaded shaft. Use a 2.5 mm key to tighten the shaft. See [Figure 360](#).

**Figure 360.** Assembling Tee Connector and Mounting Bracket (3)

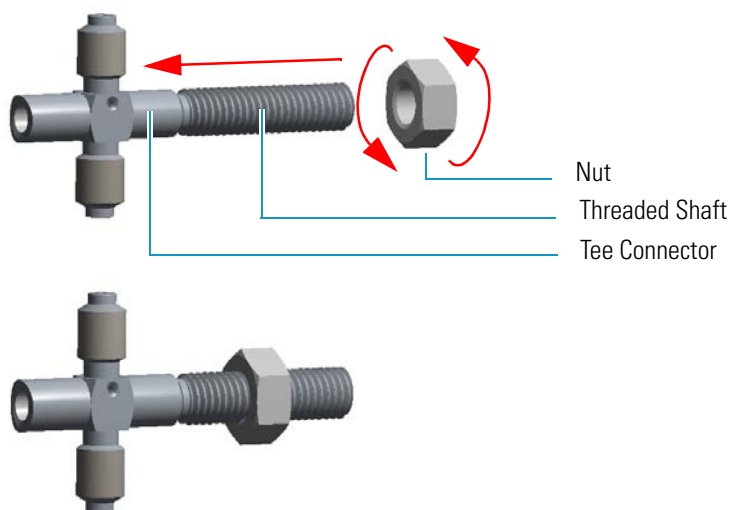

- Screw the Tee connector into clamp. Turn the threaded shaft until it protrudes the clamp, then turn the nut until it grips the clamp. See [Figure 361](#).

**Figure 361.** Assembling Tee Connector and Mounting Bracket (4)

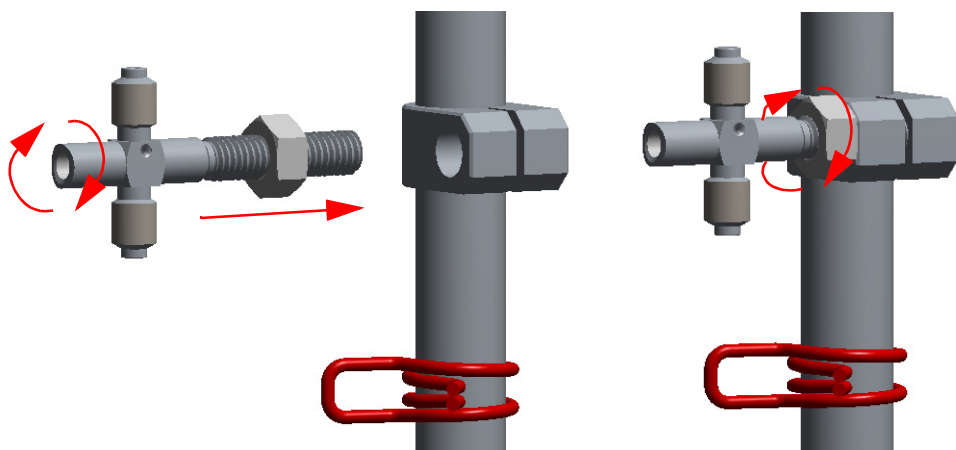

5. Move the Tee connector support on the lower part of the mounting bracket, then finger-tighten the nut up to grip the clamp without overtightening. This allow to adjust the position of the Tee connector during the installation of backflush line, precolumn, and capillary column.
6. At this point mount the mounting bracket into the oven of the GC. See [“To install the mounting bracket”](#) on [page 367](#).

❖ **To install the mounting bracket into the GC oven**

Install the mounting bracket near the front of the GC oven on the right-hand side. This will keep the mounting bracket out of the way of the column.


Proceed as follows:

1. Make sure that the nut at the top end of the mounting bracket are removed. See [Figure 358](#) on [page 361](#).
2. Line up the bottom of the mounting bracket into two holes on the bottom side of the GC oven.
3. Hold the mounting bracket directly upright, and twist the top end of it so that it is securely attached to the GC oven holes.
4. Let out the top of the mounting bracket until it is long enough to be secured into two holes on the bottom of the GC oven.
5. Replace the nut.
6. Tighten the nut until the mounting bracket is securely attached to the top and bottom of the GC oven.

[Figure 362](#) shows the result of the operation.



**Figure 362.** Mounting Bracket Installed in the GC

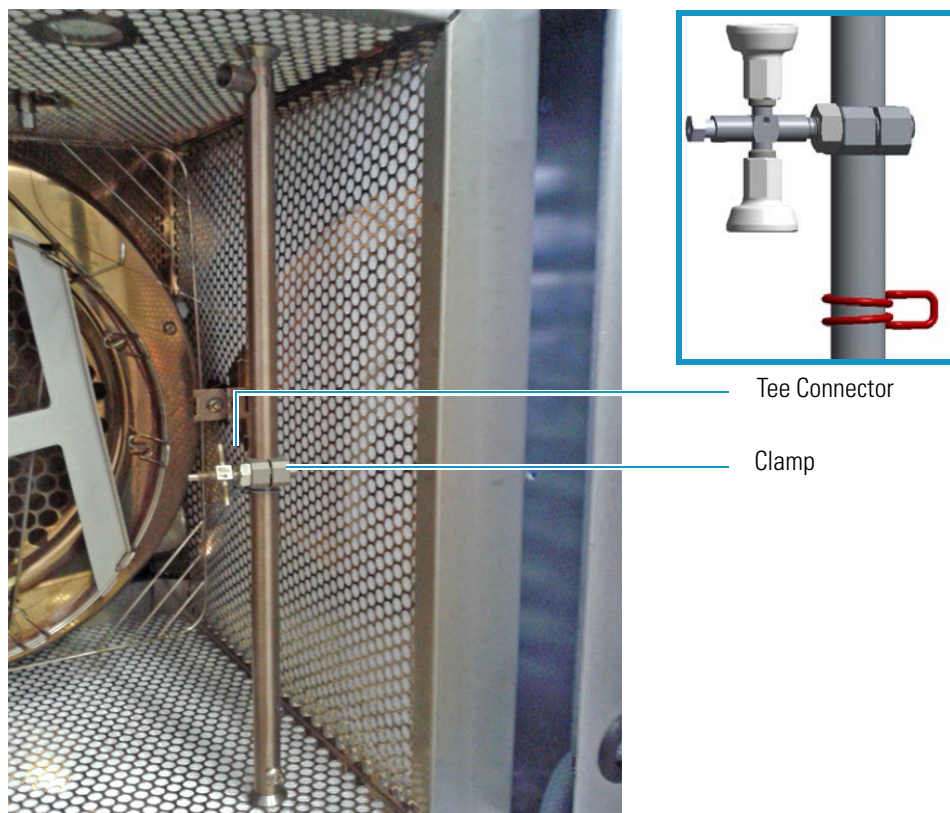


7. At this point connect backflush line, precolumn, and capillary column to the Tee connector. See [“To connect backflush line, precolumn, and capillary column”](#) on [page 369](#).

❖ **To connect backflush line, precolumn, and capillary column**

The result of this operation is schematically shown in [Figure 368](#).

**Figure 363.** Connecting Backflush Line, Precolumn, and Capillary Column

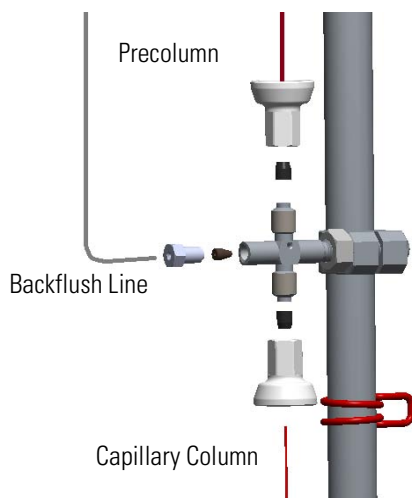

**Note** Could be necessary to adjust the position of the Tee connector when installing backflush line, precolumn, and column to stress as minimal as possible.

1. Connect the backflush line metal tube coming from the bottom of the SSLBKF injector module or PTVBKF injector module to the Tee connector.

**Tip** The backflush line must be bent as required by the position of the Tee connector on the mounting bracket.

- a. Slide the retaining nut and the ferrule onto the end of the backflush line.
  - b. Insert the backflush line into the Tee connector.
  - c. Finger-tighten the retaining nut until it starts to grip the backflush line.
  - d. Use the 5 mm wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal without overtighten.
2. Connect the precolumn to the Tee connector.


**Note** Always use original Thermo Fisher Scientific precolumns.  
The use of precolumns not meeting the specifications of our products does not ensure a good operation of the instrument and may even compromise the analytical results.

Thermo Fisher Scientific precolumn (2 m; 10 m) are mounted on the appropriate rack and identified with a relevant label fixed on the rack.

- a. Place the precolumn on the column support.
  - b. Slide the retaining nut and the 0.8 mm Vespel/graphite ferrule onto the precolumn with the flat end facing towards the Tee connector. Be careful to avoid damaging the graphite ferrule when inserting the precolumn. Always use a new ferrule of the correct diameter.
  - c. Cut 1 cm from the precolumn end.
  - d. Insert the precolumn into the Tee connector in the side with the dot that it must be normally be turned towards the ceiling of the GC oven.
  - e. Finger-tighten the precolumn retaining nut until it starts to grip the precolumn.
  - f. Use the 1/4-in. wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal without overtighten.
3. Connect the precolumn to the injector base.
    - a. Slide the retaining nut and the graphite ferrule onto the fused silica precolumn (0.53 mm ID; 2 m length) with the bevelled end facing towards the injector.
    - b. Cut 1 cm from the precolumn.
    - c. Insert the precolumn into the injector and slide the ferrule up to the injector base. Insert the precolumn about 3 cm into the bottom of the injector.

- d. Finger-tighten the precolumn retaining nut until it starts to grip the precolumn.
  - e. Use the 6 mm wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal without overtighten.
4. Connect the capillary column to the Tee connector and to the precolumn.  
See [Figure 369](#).

**Figure 364.** Precolumn and Capillary Column Connection

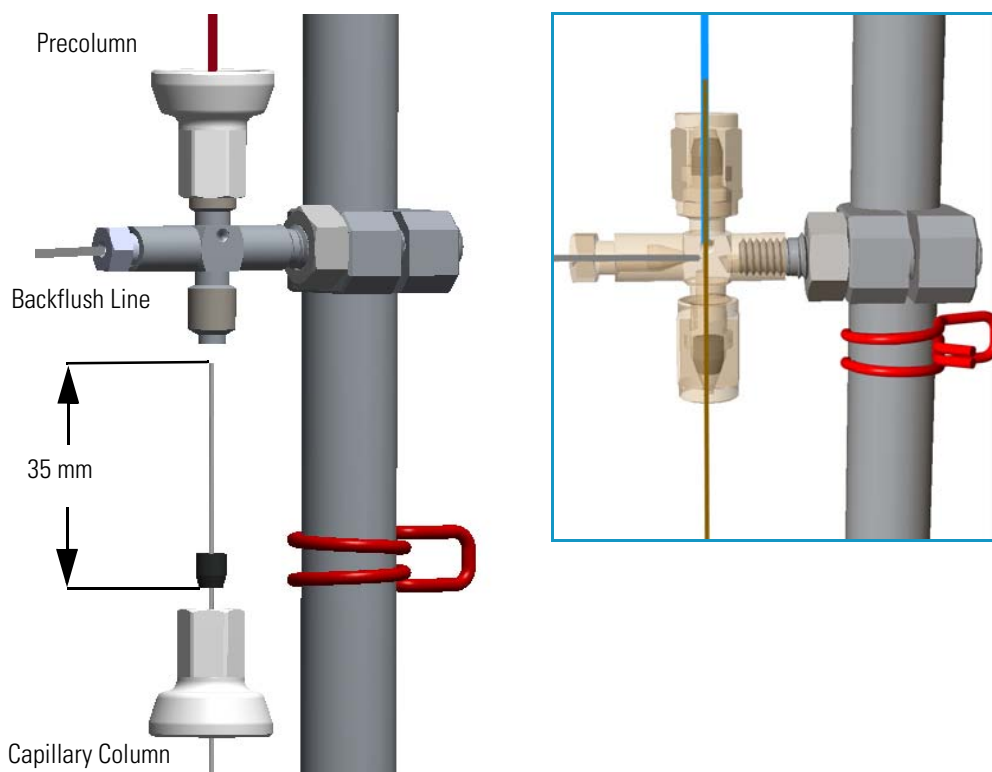

- a. Place the column on the column support.
- b. Slide the proper retaining nut and Vespel/graphite ferrule onto the column with the flat end facing towards the Tee connector. Be careful to avoid damaging the ferrule when inserting the column.
- c. Cut 1 cm from the column end.
- d. Use typewriter correction fluid or a felt-tipped pen to mark the position of the ferrule bottom 35 mm from the end of the column.
- e. Carefully insert the column into the precolumn up to reach the marked position. The capillary column end must protrudes the Tee connector; if not, repeat the operations.
- f. Finger-tighten the column retaining nut until it starts to grip the column.
- g. Use the 1/4-in. wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal.


5. Connect the analytical column to the detector base.
  - a. Slide the retaining nut and the graphite ferrule onto the capillary column with the bevelled end facing the detector base.

Be careful to avoid damaging the graphite ferrule when inserting the column.
  - b. Cut 2–3 cm from the column end.
  - c. Insert the column into the detector base body and slide the ferrule up to the detector base.
  - d. Finger-tighten the column retaining nut until it starts to grip the column.
  - e. Push the column through the detector base according to the detector in use.
  - f. Use the 6 mm wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal without overtighten.

❖ **To install of a second backflush system**

If a second SSL/PTV Backflush system is installed, you can connect both the systems by using two Tee connectors mounted on the same mounting bracket as schematically shown in [Figure 370](#).

**Figure 365.** Schematically Connection of a Second Backflush System

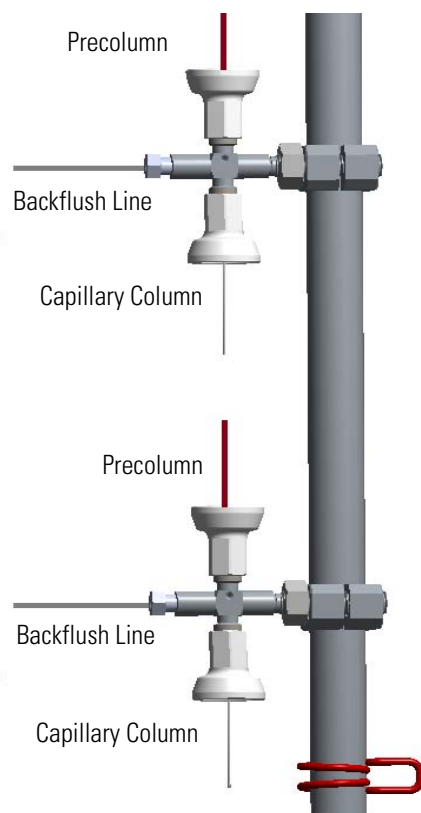

**Tip** Turn the Tee connector as required to allow the installation of the two precolumns.

## Connecting a SSL/PTV Backflush System for High Temperature

This section provides instructions for connecting your SSLBKF/PTVBKF injector modules with precolumn and analytical capillary column into the oven for application at high temperature. For the purpose, the kit PN 19050254 is required. Refer to the chapters SSLBKF Injector Module and PTVBKF Injector Module in the *TRACE 1300/TRACE 1310 User Guide*. This system operates at a temperature more than 300 °C.

The operation consists of the following operating procedures:

- “To install the mounting bracket” on page 367
- “To connect backflush line, precolumn, and capillary column” on page 369
- “To install of a second backflush system” on page 373



**CAUTION** Before starting make sure the SSLBKF/PTVBKF injector module is correctly installed into its seat.

### ❖ To install the mounting bracket

Install the mounting bracket near the front of the GC oven on the right-hand side. This will keep the mounting bracket out of the way of the column.

Proceed as follows:

1. Remove the nut at the top end of the mounting bracket. This nut secures the two sections of the mounting bracket. See [Figure 366](#).

**Figure 366.** Mounting Bracket Parts




2. Push the top of the mounting bracket down.
3. Line up the bottom of the mounting bracket into two holes on the bottom side of the GC oven.
4. Hold the mounting bracket directly upright, and twist the top end of it so that it is securely attached to the GC oven holes.
5. Let out the top of the mounting bracket until it is long enough to be secured into two holes on the bottom of the GC oven.
6. Replace the nut.

## 8 Installing Optional Kits

Connecting a SSL/PTV Backflush System for High Temperature

7. Tighten the nut until the mounting bracket is securely attached to the top and bottom of the GC oven.


**Figure 367.** Mounting Bracket Installed in the GC

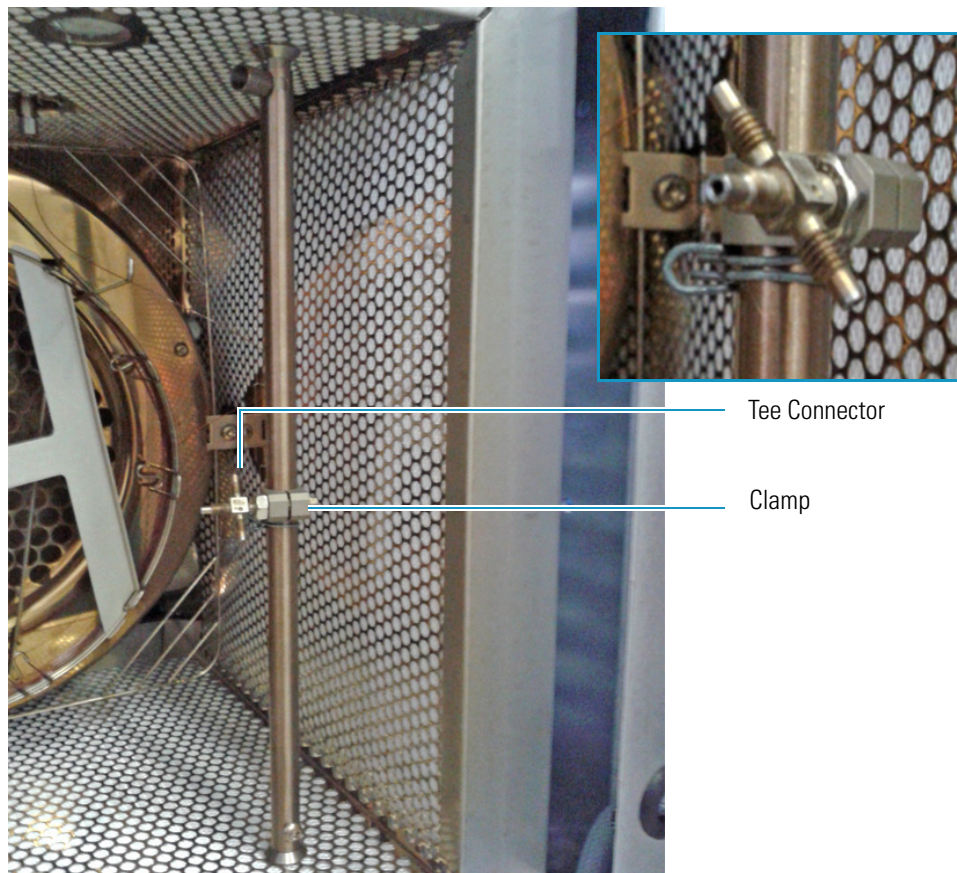




**Note** If the Tee connector should not be assembled into its clamp, proceed as follows:

1. Completely screw the nut on the threaded shaft of the Tee connector.

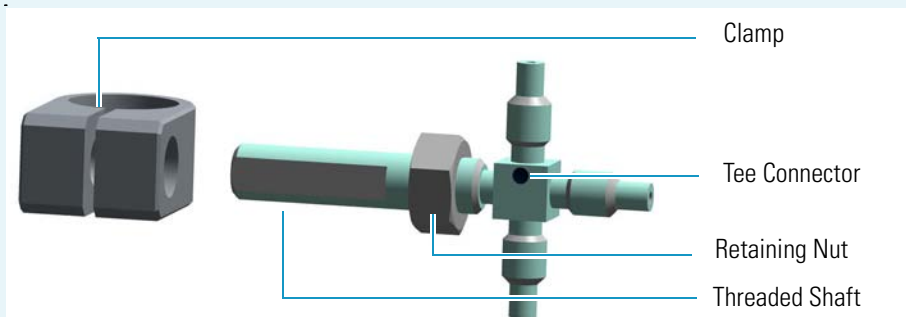


Note that the clamp has two holes: the smaller one is threaded while the bigger other is not threaded.

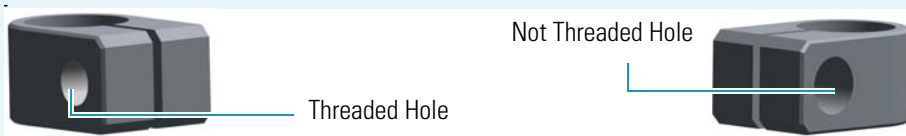

2. Pay attention to insert the threaded shaft of the Tee connector into the clamp through the bigger hole not threaded up to grip the threaded hole.



3. Turn the Tee connector until the threaded shaft begins to protrude the clamp.

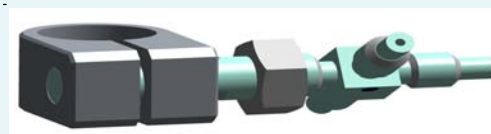


4. Move the Tee connector support on the lower part of the mounting bracket, then finger-tighten the nut up to grip the clamp without overtightening. This allow to adjust the position of the Tee connector during the installation of backflush line, precolumn, and capillary column.

❖ **To connect backflush line, precolumn, and capillary column**

The result of this operation is schematically shown in [Figure 368](#).

**Figure 368.** Connecting Backflush Line, Precolumn, and Capillary Column

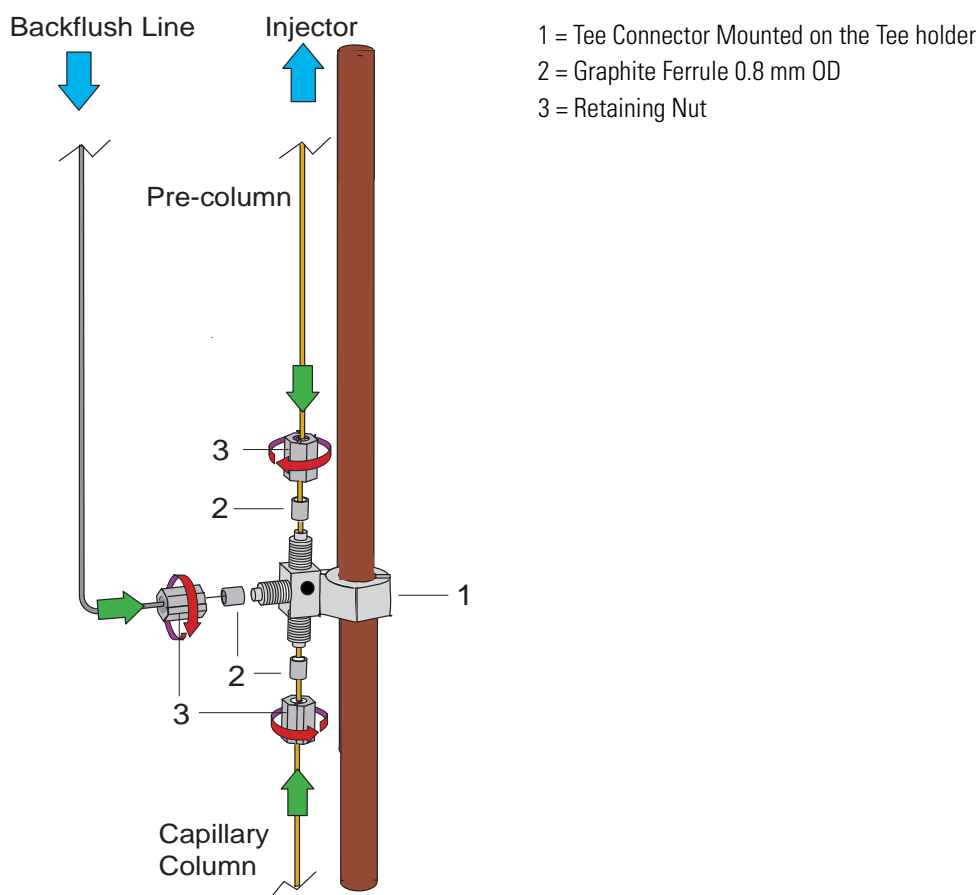

**Note** Could be necessary to adjust the position of the Tee connector when installing backflush line, precolumn, and column to stress as minimal as possible.

1. Connect the backflush line metal tube coming from the bottom of the SSLBKF injector module or PTVBKF injector module to the Tee connector.

**Tip** The backflush line must be bent as required by the position of the Tee connector on the mounting bracket.

- a. Slide the retaining nut and the 1.0 mm graphite ferrule onto the end of the backflush line with the bevelled end facing towards the Tee connector. Be careful to avoid damaging the graphite ferrule when inserting the tube.
- b. Shortly activate the backflush with carrier open to eliminate possible graphite pieces before connecting the Tee connector.
- c. Insert the backflush line into the Tee connector.
- d. Finger-tighten the retaining nut until it starts to grip the backflush line.



- e. Use the 6 mm wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal without overtighten.
2. Connect the precolumn to the Tee connector.


**Note** Always use original Thermo Fisher Scientific precolumns.

The use of precolumns not meeting the specifications of our products does not ensure a good operation of the instrument and may even compromise the analytical results.

Thermo Fisher Scientific precolumn (2 m; 10 m) are mounted on the appropriate rack and identified with a relevant label fixed on the rack.

- a. Place the precolumn on the column support.
  - b. Slide the retaining nut and the 0.8 mm graphite ferrule onto the precolumn with the bevelled end facing towards the Tee connector. Be careful to avoid damaging the graphite ferrule when inserting the precolumn. Always use a new ferrule of the correct diameter.
  - c. Cut 1 cm from the precolumn end.
  - d. Insert the precolumn into the Tee connector in the side with the dot that it must be normally be turned towards the ceiling of the GC oven.
  - e. Finger-tighten the precolumn retaining nut until it starts to grip the precolumn.
  - f. Use the 6 mm wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal without overtighten.
3. Connect the precolumn to the injector base.
  - a. Slide the retaining nut and the graphite ferrule onto the fused silica precolumn (0.53 mm ID; 2 m length) with the bevelled end facing towards the injector.
  - b. Cut 1 cm from the precolumn.
  - c. Insert the precolumn into the injector and slide the ferrule up to the injector base. Insert the precolumn about 3 cm into the bottom of the injector.
  - d. Finger-tighten the precolumn retaining nut until it starts to grip the precolumn.
  - e. Use the 6 mm wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal without overtighten.
4. Connect the capillary column to the Tee connector and to the precolumn.  
See [Figure 369](#).

**Figure 369.** Precolumn and Capillary Column Connection

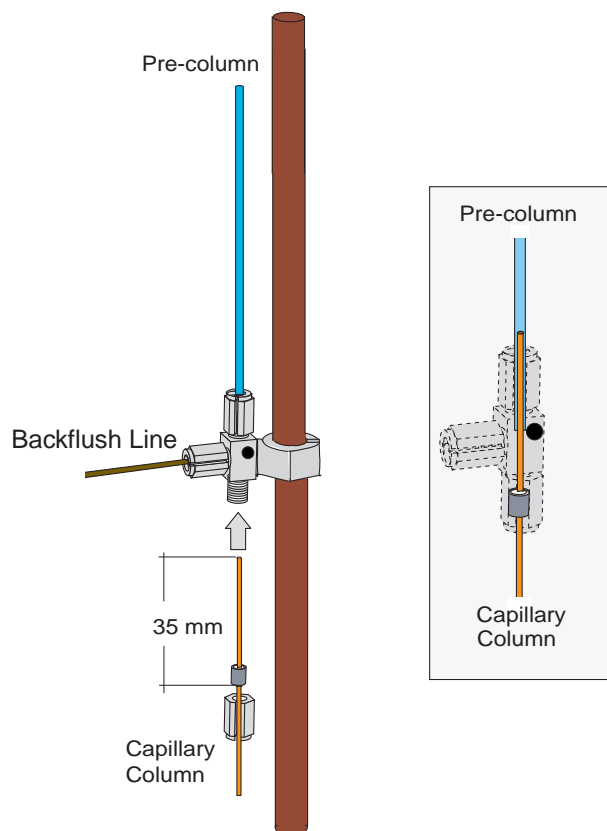

- a. Place the column on the column support.
  - b. Slide the proper retaining nut and graphite ferrule onto the column with the bevelled end facing towards the Tee connector. Be careful to avoid damaging the graphite ferrule when inserting the column.
  - c. Cut 1 cm from the column end.
  - d. Use typewriter correction fluid or a felt-tipped pen to mark the position of the ferrule bottom 35 mm from the end of the column.
  - e. Carefully insert the column into the precolumn up to reach the marked position. The capillary column end must protrudes the Tee connector; if not, repeat the operations.
  - f. Finger-tighten the column retaining nut until it starts to grip the column.
  - g. Use the 6 mm wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal.
5. Connect the analytical column to the detector base.
    - a. Slide the retaining nit and the graphite ferrule onto the capillary column with the bevelled end facing the detector base.


Be careful to avoid damaging the graphite ferrule when inserting the column.

- b. Cut 2–3 cm from the column end.
- c. Insert the column into the detector base body and slide the ferrule up to the detector base.
- d. Finger-tighten the column retaining nut until it starts to grip the column.
- e. Push the column through the detector base according to the detector in use.
- f. Use the 6 mm wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal without overtighten.

❖ **To install of a second backflush system**

If a second SSL/PTV Backflush system is installed, you can connect both the systems by using two Tee connectors mounted on the same mounting bracket as schematically shown in [Figure 370](#).

**Figure 370.** Schematically Connection of a Second Backflush System

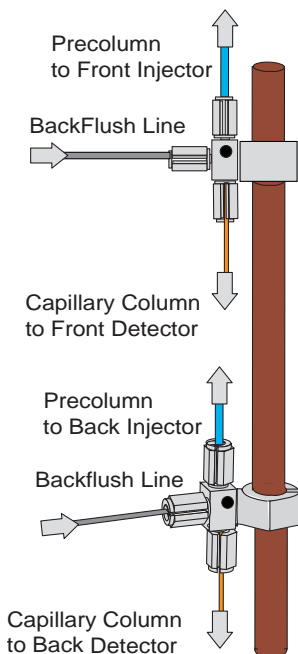

**Tip** Turn the Tee connector as required to allow the installation of the two precolumns.

## Connecting a GSV Backflush System

This sections provides instructions for connecting your Gas Sampling Valve module with precolumn and analytical capillary column into the oven of the TRACE 1300/1310 GC. Refer to the chapters Gas Sampling Valve (GSV) Module in the *TRACE 1300/TRACE 1310 User Guide*.

The Gas Sampling Valve backflush system is provided with the kit PN 19050764.


❖ **To connect backflush line, precolumn, and analytical capillary column**



**CAUTION** Before starting make sure the Gas Sampling Valve (GSV) module is correctly installed into its seat.

The result of this operation is schematically shown in the example of [Figure 371](#).

**Figure 371.** Connecting GSV Backflush Line, Precolumn, and Analytical Capillary Column

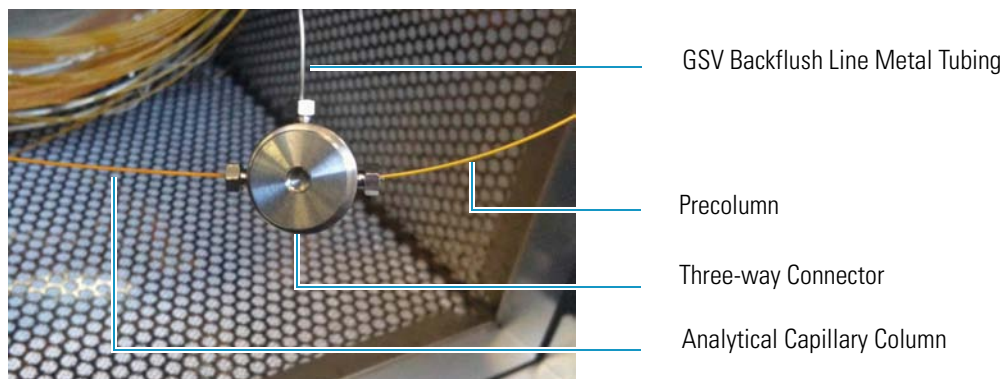

**Note** Could be necessary to adjust the position of the three-way connector when installing backflush line, precolumn, and capillary column to stress as minimal as possible.

When the GSV module is inserted into its position on the upper deck of the GC, the backflush line metal tubing protrudes into the GC oven.

1. Open the front door of the GC.
2. Remove nut and ferrule from the metal tubing of the GSV module.
3. Connect the backflush line metal tubing coming from the bottom of the Gas Sampling Valve module to the three-way connector.
  - a. Slide the 1/32-in. nut and the Vespel/Graphite ferrule PN 29003428 onto the end of the backflush line with the bevelled end facing towards the three-way connector.
  - b. Insert the backflush line into the central position of the three-way connector.
  - c. Slide the retaining nut onto the backflush line through its side cut.
  - d. Use the 5 mm wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal without overtighten.
4. Connect the precolumn to the Tee connector.
  - a. Place the precolumn on the column support.
  - b. Slide the 1/32-in. nut and the proper Vespel/Graphite ferrule onto the precolumn with the bevelled end facing towards the three-way connector. Always use a new ferrule of the correct diameter.
  - c. Cut 1 cm from the precolumn end.

- d. Insert the precolumn in one of the two external positions of the three-way connector.
  - e. Finger-tighten the precolumn retaining nut until it starts to grip the precolumn.
  - f. Use the 5 mm wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal without overtighten.
5. Connect the precolumn to the injector base.
- a. Slide the M6 nut and the graphite ferrule onto the fused silica precolumn with the bevelled end facing towards the injector.
  - b. Cut 1 cm from the precolumn.
  - c. Insert the precolumn into the injector and slide the ferrule up to the injector base. Insert the precolumn till it reaches the mechanical stop, than drawn it back of 1 mm.
  - d. Finger-tighten the precolumn retaining nut until it starts to grip the precolumn.
  - e. Use the 1/4-in. wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal without overtighten.
6. Connect the capillary column to the three-way connector.
- a. Place the capillary column on the column support.
  - b. Slide the 1/32-in. nut and the proper Vespel/Graphite ferrule onto the capillary column with the bevelled end facing towards the three-way connector.
  - c. Cut 1 cm from the end of the capillary column.
  - d. Insert the end of the capillary column in the free position of the three-way connector.
  - e. Finger-tighten the column retaining nut until it starts to grip the column.
  - f. Use the 5 mm wrench to tighten the retaining nut. Use enough pressure necessary to obtain a good seal.
7. Connect the capillary column to the detector base.
- a. Slide the nut and the graphite ferrule onto the capillary column with the bevelled end facing the detector base.
  - b. Be careful to avoid damaging the graphite ferrule when inserting the capillary column.
  - c. Cut 2–3 cm from the end of the capillary column.
  - d. Insert the capillary column into the detector base body and slide the ferrule up to the detector base.
  - e. Finger-tighten the capillary column retaining nut until it starts to grip the column.
  - f. Push the capillary column through the detector base according to the detector in use.
  - g. Use the 1/4-in. wrench to tighten the retaining nut.
  - h. Use enough pressure necessary to obtain a good seal without overtighten.

## Installing the NoVent Microfluidics

This section provides instructions for installing the NoVent™ Microfluidics on your TRACE 1300 or TRACE 1310 GC, and your ISQ or TSQ 8000 Series mass spectrometer.


### Contents

- “Connecting the NoVent Microfluidics Module to the TRACE 1300/1310” on page 376
- “Installing the Mounting Bracket” on page 379
- “Preparing the NoVent Microfluidics Restrictor Tubing” on page 381
- “Attaching the Ferrule and Nut to the GC Column” on page 384
- “Attaching the New Tubing to the Transfer Line” on page 385
- “Connecting the Capillaries to the Microfluidics Splitter” on page 390
- “Configuring the Post-Column” on page 391
- “Using the Module” on page 394

## Connecting the NoVent Microfluidics Module to the TRACE 1300/1310

According to the configuration of your TRACE 1300 or TRACE 1310 GC, a dummy module is present in the free site where the detector module is not installed. If no dummy module is present, and no vacant module position is present on either the TRACE 1300 or TRACE 1310 GC or on the optional gas valve oven (if you have one), you must take one of the existing and installed modules off in order to install the NoVent Microfluidics.

**Figure 372.** Location of Detector Modules on your TRACE 1300/TRACE 1310 GC

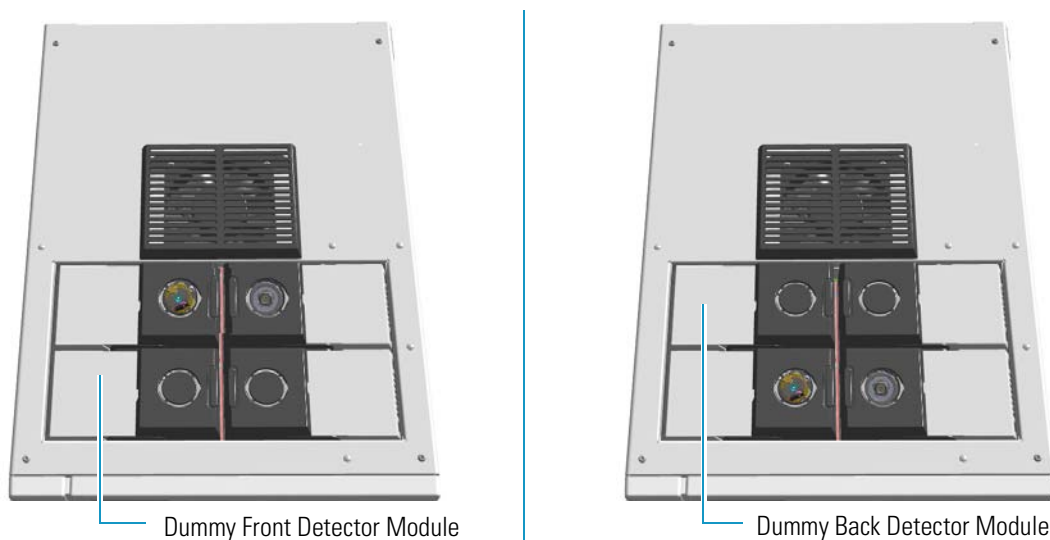



**ATTENTION** Where a dummy module is installed, the gas connection is blocked by a plug.

❖ **To connect the NoVent Microfluidics module to the TRACE 1300/1310 GC**

1. Cool the oven, injector or injectors, transfer line, ion source, and any installed GC detectors to room temperature and shut down the GC.




**WARNING - BURN HAZARD:** The injectors, detectors, oven, and transfer line may be hot. Allow them to cool to room temperature before touching them.

2. Push down on the power switch to power off the GC.
3. Shut down the TSQ 8000 Series or ISQ Series mass spectrometer using the software. See your mass spectrometer's user guide for more information. The heaters and the turbo-molecular pump power off.
4. If you are using hydrogen as a carrier gas, unscrew the hydrogen safety screw on the front door of the mass spectrometer.




**WARNING - FIRE HAZARD:** If you are using hydrogen, do NOT reach over the top of the instrument to power it off. Instead, reach around the right side or go to the back of the instrument and flip down the power switch.

**Figure 373.** Finding the Vent Valve on the Mass Spectrometer

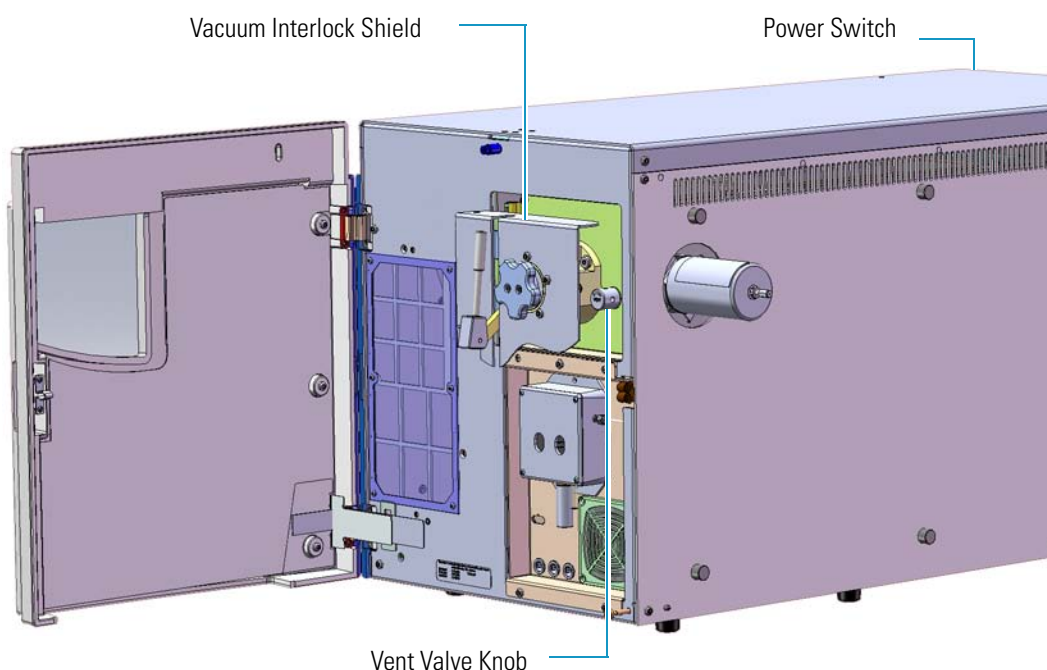

5. Open the front door of the mass spectrometer.
6. Look behind the right side of the vacuum interlock shield and twist the vent valve knob one and a half times in a counter-clockwise direction to open the vent valve.
7. Wait five minutes for the mass spectrometer to vent.
8. Open the front door of the GC and loosen the transfer line nut. Then pull the column back (into the oven) about 5 cm to ensure the column is no longer in the ion source.
9. Close all GC gas supplies. Shut the carrier gas supply off at its source, such as the tank.
10. Remove the dummy module from the position where the NoVent Microfluidics will be installed. The most convenient location for accessing the tubing connections is the front position.
  - a. Open the module flap cover.
  - b. Using a T20 Torxhead screwdriver, unscrew and remove the two captive fixing screws.
  - c. Keeping the dummy module flap cover open, lift up the module from its seat in the detector housing. Place the dummy module on a clean surface.
  - d. Remove the gas plug by unscrewing its fixing screw using a T20 Torxhead screwdriver.

**Figure 374.** Removing the Detector Gas Plug

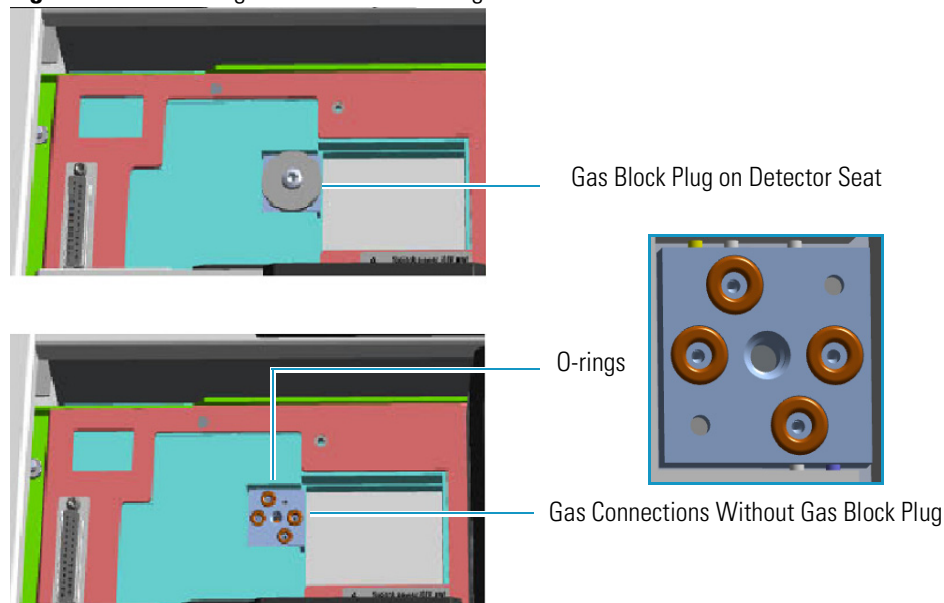

**WARNING** Make sure all four o-rings are correctly seated on the gas connection. Do not install the module if the o-rings are missing.



11. Install the No-Vent module.
  - a. Open the module flap cover.
  - b. Keeping the module flap cover open, place the module in its seat. Be sure to insert the 25-pin male connector, on the bottom of the module, into the 25-pin female connector on the detector seat of the detector housing.
  - c. Use a T20 screwdriver to tighten the three captive fixing screws without overtightening.

**Note** Tighten the center screw first and then secure the side screws.

- d. Close the module flap cover.
12. Power on the GC.
13. Switch the module on so that the LED light glows green and start the carrier gas flowing to purge the line.

**Note** This module requires a constant carrier gas pressure of 60 psig (410 kPa).

## Installing the Mounting Bracket


Install the mounting bracket near the front of the GC oven on the right-hand side. This will keep the mounting bracket out of the way of the column.

**Note** If you already have a mounting bracket installed in your GC with either the SSL backflush or PTV backflush system, skip this step.

### ❖ To install the mounting bracket

1. Remove the nut at the top end of the mounting bracket. This nut secures the two sections of the mounting bracket. See [Figure 375](#).


**Figure 375.** Mounting Bracket Parts

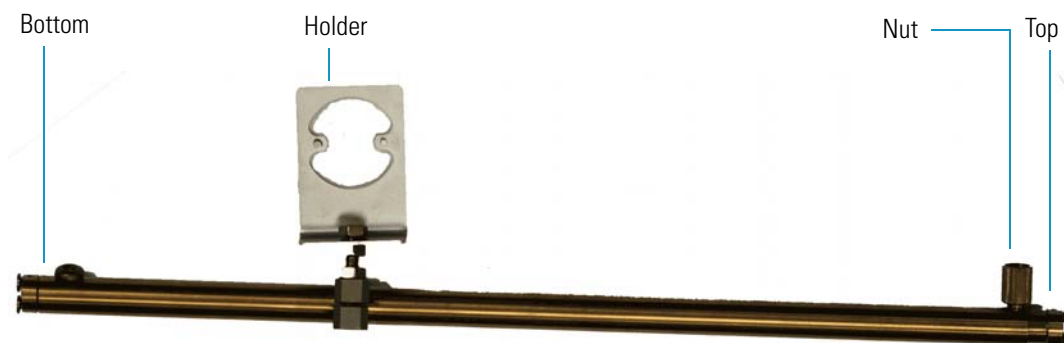

2. Push the top of the mounting bracket down.
3. Line up the bottom of the mounting bracket into two holes on the bottom side of the GC oven.


4. Hold the mounting bracket directly upright, and twist the top end of it so that it is securely attached to the GC oven holes.
5. Let out the top of the mounting bracket until it is long enough to be secured into two holes on the bottom of the GC oven.
6. Loosely replace the nut, then rotate the top and bottom in opposite directions.
7. Tighten the nut until the mounting bracket is securely attached to the top and bottom of the GC oven.

**Figure 376.** Mounting Bracket Installed in the GC

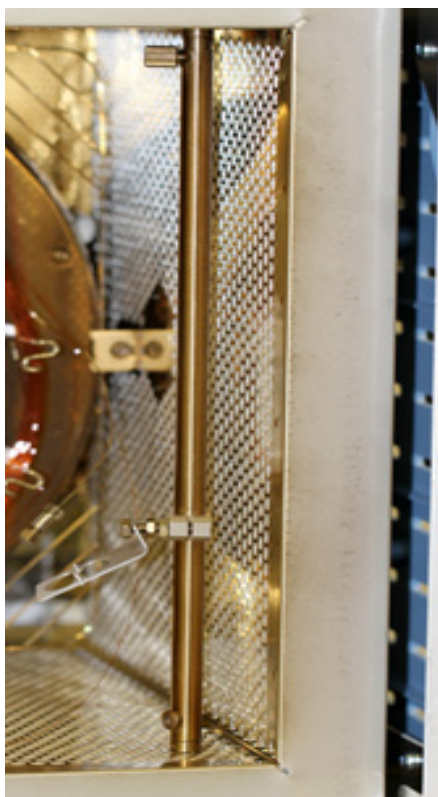

8. Attach the Thermo Scientific™ microfluidics splitter to the mounting bracket. The microfluidics splitter snaps into place. The small hole should be positioned at the top. See [Figure 377](#) for the correct orientation.

**Figure 377.** Correct Orientation for the Microfluidics Splitter

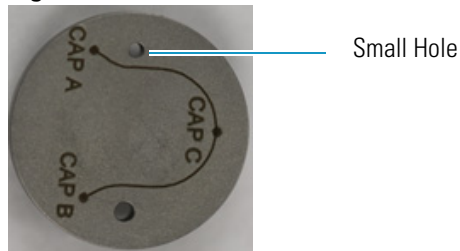

9. Position the holder for the microfluidics splitter upright and use two 7 mm wrenches to secure it. See [Figure 378](#).

**Figure 378.** Securing the Holder for the Microfluidics Splitter



## Preparing the NoVent Microfluidics Restrictor Tubing

If you are using helium as a carrier gas, use the 30 cm length of 75  $\mu$ m fused silica tubing in your kit (two are supplied) to connect the NoVent Microfluidics module to the splitter. If you are using hydrogen as a carrier gas, use the 80 cm length of 75  $\mu$ m fused silica tubing in your kit (one is supplied) to connect the NoVent Microfluidics module to the splitter. This tubing acts as a flow restrictor to limit the helium flow from the NoVent microfluidics module. You must swage a SilTite™ ferrule to this tubing. You will need the following materials to connect the tubing to new ferrule and nut. They are all included in the complete kit.


- SilTite FingerTite™ pre-swage tool—to pre-swage the SilTite FingerTite ferrule to the GC module
- SilFlow™ pre-swage tool and SilFlow FingerTite tool—to pre-swage the SilFlow ferrule to the microfluidics splitter
- SilFlow nut—to connect the tubing to the microfluidics splitter
- SilTite ferrule (1)
- SilFlow ferrules (3)

**Note** If you are using a 0.25 mm id column, use a 0.4 mm SilFlow ferrule and a 0.4 mm SilFlow FingerTite jig. If you are using a 0.32 mm id column, use a 0.5 mm SilFlow ferrule and a 0.5 mm SilFlow FingerTite jig.

1. Use a scoring wafer to cut the end off the tubing.
2. Use the SilTite FingerTite pre-swage tool in your kit to pre-swage the SilTite ferrule to your tubing. See [Figure 379](#). The tubing should extend just past the tool when pre-swaging the ferrule. After the ferrule is swaged, remove the tool and confirm the ferrule does not slide on the tubing.

Then cut the tubing to about 2 mm past the tip of the ferrule. See [Figure 379](#).

**Figure 379.** Pre-swaging the Ferrule to the 30 cm Fused Silica Tubing

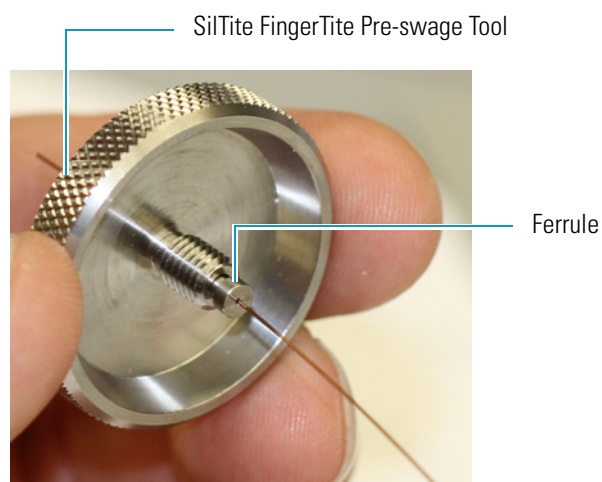

3. Insert the ferrule and the tubing into the NoVent microfluidics module. See [Figure 380](#).

**Figure 380.** Securing the 30 cm Fused Silica Restrictor to the No-Vent Microfluidics Module



4. Add the FingerTite screw to the module and tighten.
5. Now you must prepare the other end of the tubing to connect to the microfluidics 3-port splitter. It is much easier to attach the ferrule to the tubing while it is outside the GC. You will connect it to the splitter later.
6. Position a SilFlow nut and SilFlow ferrule onto the tubing as shown in [Figure 381](#).

**Figure 381.** Positioning the SilFlow Nut and SilFlow Ferrule Correctly on the Tubing



7. Use a scoring wafer to cut the tubing after inserting it through the ferrule. See [Figure 382](#). Then use the SilFlow pre-swage tool to secure the ferrule into position. When done properly, the tubing will extend slightly past the tip of the ferrule.


It is important to use the SilFlow pre-swage tool in order to prevent crushing the tip of the fused silica. See [Figure 382](#).

**Figure 382.** Cutting the Tubing with the Scoring Wafer



8. Place the column and ferrule into SilFlow pre-swage tool until the tubing reaches the bottom of the tool. When done properly, the tubing will extend slightly past the tip of the ferrule. It is important to use the pre-swage tool in order to prevent crushing the tip of the fused silica. See [Figure 383](#).

**Figure 383.** Inserting the Column and Ferrule into the SilFlow Pre-Swage Tool

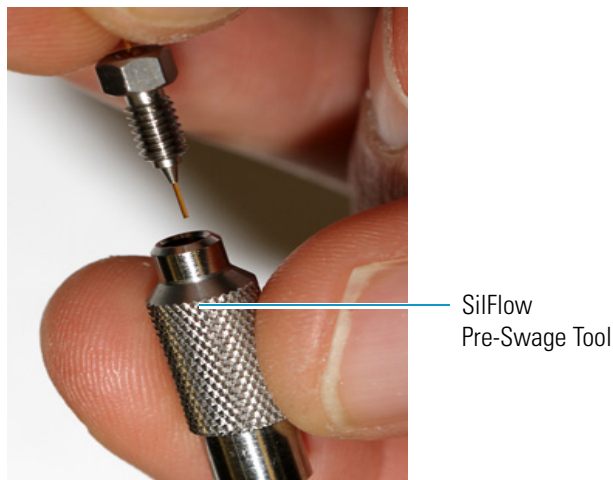

9. Use the SilFlow FingerTite tool to swage the ferrule to the tubing. Be sure to keep the tip of the fused silica bottomed out in the pre-swage tool. See [Figure 384](#).

**Figure 384.** Swaging the Ferrule Using the SilFlow FingerTite Tool



10. Remove the jig, and lay the tubing carefully on the bottom of the GC until you are ready to connect it to the microfluidics splitter.

## Attaching the Ferrule and Nut to the GC Column

You will need the following materials to connect the column to new ferrule and nut. They are all included in the complete kit.


- SilFlow FingerTite tool
- SilFlow pre-swage tool
- SilFlow nut
- SilFlow ferrules

**Note** If you are using a 0.25 mm id column, use a 0.4 mm SilFlow ferrule and a 0.4 mm SilFlow pre-swage tool. If you are using a 0.32 mm id column, use a 0.5 mm SilFlow ferrule and a 0.5 mm SilFlow pre-swage tool.

### ❖ To add the new ferrule and nut to the column

1. Attach the ferrule to the column while it is outside the GC. You will connect it to the module later.
2. Position the SilFlow nut and SilFlow ferrule onto the column as shown in [Figure 385](#).

**Figure 385.** Positioning the Nut and Ferrule Correctly on the Column

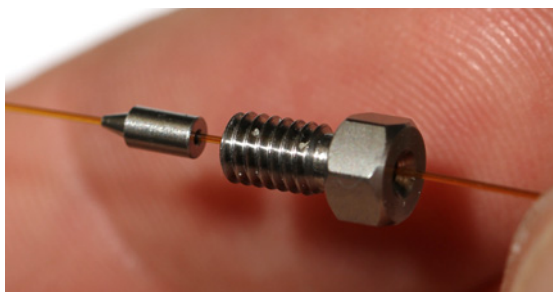

3. Use a scoring wafer to cut the column after inserting it through the ferrule. Then use the appropriate pre-swage tool to secure the ferrule into position. When done properly, the tubing will extend slightly past the tip of the ferrule. It is important to use the pre-swage tool in order to prevent crushing the tip of the fused silica. See [Figure 386](#).

**Figure 386.** Cutting the Column with the Scoring Wafer





4. Place the column and ferrule into the SilFlow pre-swage tool until the column reaches the bottom of the jig. See [Figure 387](#).

**Figure 387.** Inserting the Column and Ferrule into the SilFlow Pre-Swage Tool



5. Use the FingerTite tool to swage the ferrule to the column. Be sure to keep the tip of the fused silica bottomed out in the pre-swage tool. See [Figure 388](#).

**Figure 388.** Swaging the Ferrule Using the SilFlow FingerTite Tool



6. Remove the pre-swage tool, and lay the column carefully on the bottom of the GC until you are ready to connect it to the microfluidics splitter.

## Attaching the New Tubing to the Transfer Line

Your kit comes with a 0.6 m and 1.2 m deactivated segments of fused silica tubing having an internal diameter of 0.17 mm. This tubing has fused ends. You must cut these ends off the tubing and attach a SilFlow ferrule before connecting it. Use this tubing to connect the microfluidics splitter to the transfer line.

**Note** Use the 60 cm tubing if you are using helium as a carrier gas. Use the 120 cm length of tubing if you are using hydrogen as a carrier gas.

When connecting the column to the transfer line, you may use either the regular transfer line nut or the spring loaded transfer line nut with the graphite Vespel<sup>®</sup> ferrule.

### ❖ To connect the column using the regular transfer line nut

1. Lower the oven temperature and allow it to cool.

2. Confirm that the MS is vented and remove the current transfer line nut and ferrule.
3. Unwind about one turn of the column from the column outlet end.


**Note** Wear clean, lint- and powder-free gloves when you handle the column and transfer line ferrule.

4. Wipe approximately 300 mm (12 in.) of the column with a tissue soaked in methanol.
5. Choose an appropriate ferrule for the outer diameter of your column.

**Note** If the maximum oven temperature in your method is  $\geq 290^{\circ}\text{C}$  ( $554^{\circ}\text{F}$ ), Thermo Fisher Scientific recommends using a spring loaded transfer line nut with a graphite Vespel® ferrule or a SilTite™ nut and ferrule. By cycling the oven at and above this temperature, expansion and contraction of the graphite Vespel® material can cause leaks in the transfer line.

6. Insert the column through the transfer line nut and ferrule, entering through the tapered end of the ferrule. Wipe the column again with a tissue soaked in methanol.

**Figure 389.** Transfer Line Nut and SilTite Ferrule Orientation

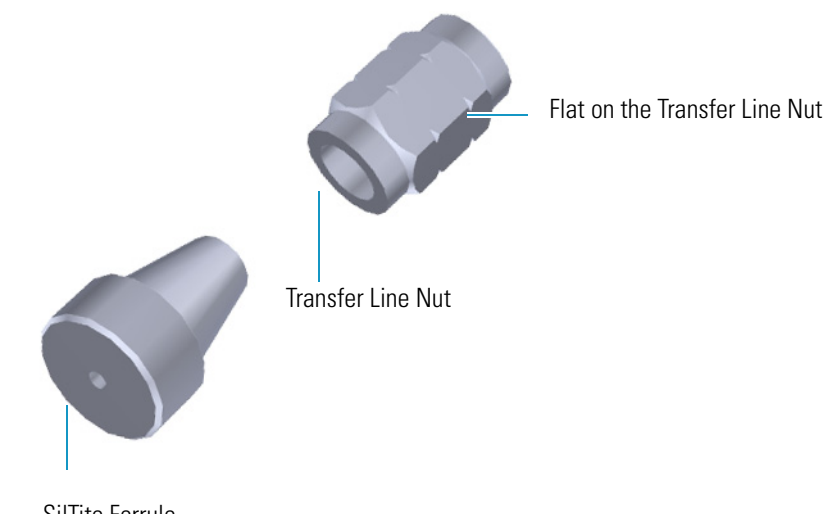

7. Insert the column into the column measuring tool (see [Figure 390](#)), which is in the ISQ Toolkit, so that it is even with the lines at the end of the column. [Figure 391](#) indicates proper positioning of the column in the tool for accurate measuring.
8. Use a scoring wafer to score and break the column. Use a magnifying glass to check for an even, flat cut. Repeat if necessary.
9. Use a 5/16 in. wrench to hold the column measuring tool steady.


**Figure 390.** Column Measuring Tool






10. While holding the column measuring tool steady, tighten the transfer line nut with a 1/4 in. wrench until the column just stops moving in the ferrule.
11. Turn the transfer line nut 1 flat backward so the column is able to move in the ferrule with slight resistance.
12. Line up the outlet of the column with the arrows on the end of the column measuring tool.

**Figure 391.** Lining Up the Column in the Column Measuring Tool

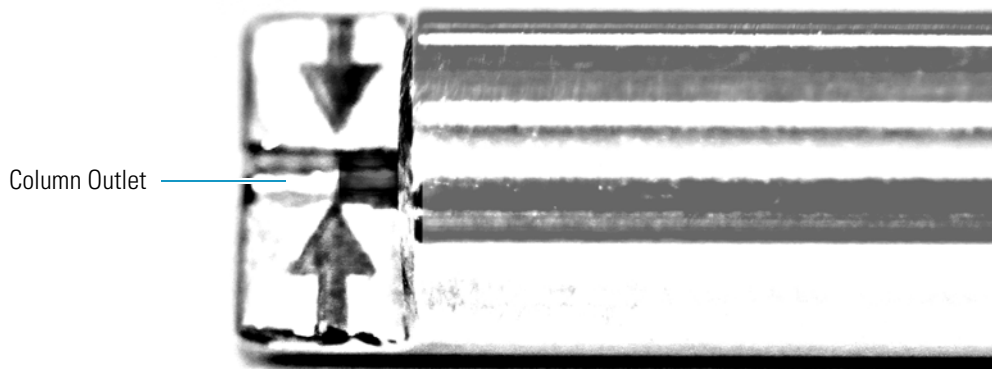


13. Place a septum with a notch cut into it behind the transfer line nut. The septum marks the place on the column where it should exit the nut.

**Figure 392.** Positioning the Septum

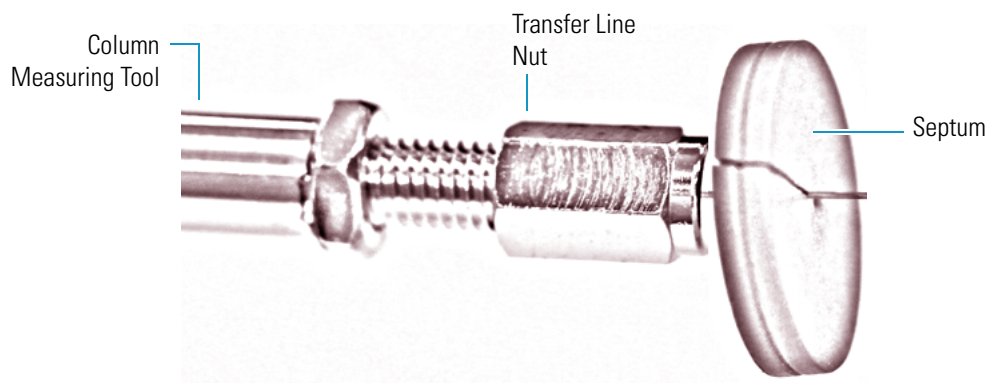


14. Pull the column back from the transfer line nut. Do not move the septum from its position on the column.

**Figure 393.** Pulling the Column Back from the Transfer Line Nut

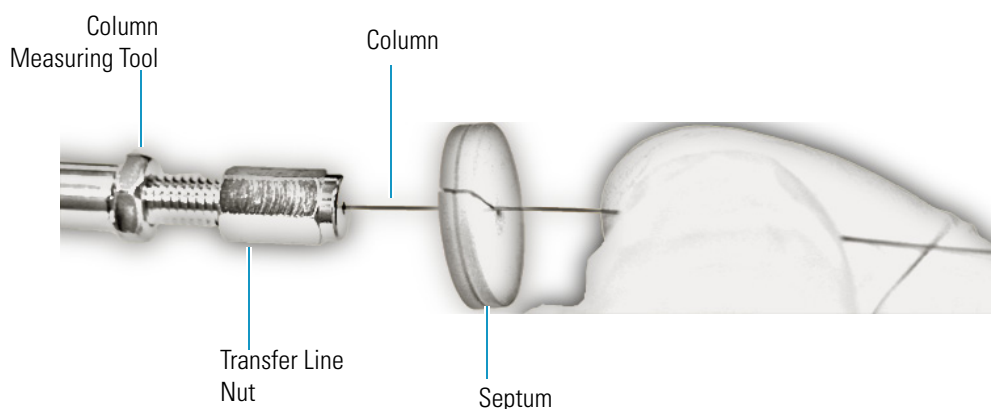


15. Loosen the transfer line nut from the column measuring tool.
16. Remove the column, transfer line nut and ferrule from the column measuring tool, making sure not to move the septum from its location on the column.
17. Insert the column into the transfer line.

**Figure 394.** Inserting the Column into the Transfer Line

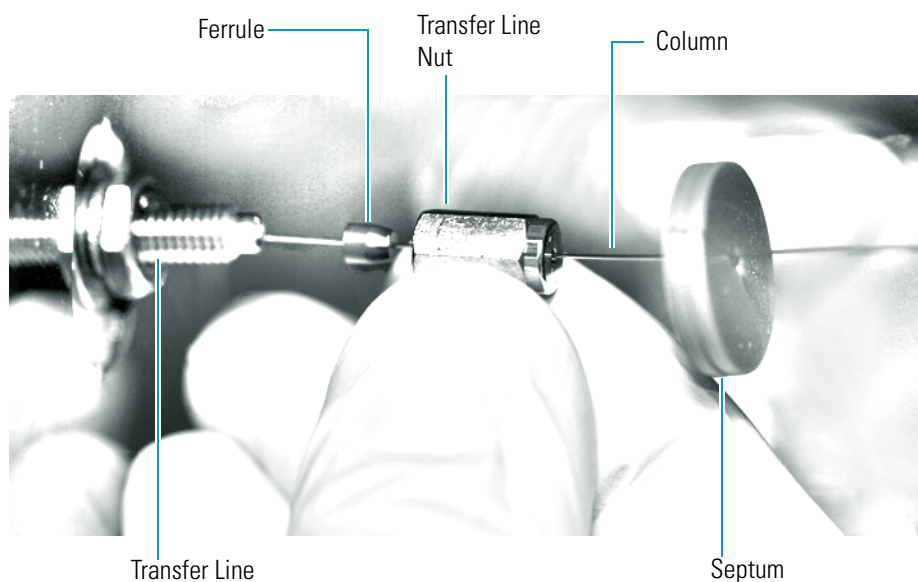

18. Tighten the transfer line nut until it is just secure enough so that you cannot move it.
19. Loosen the nut by turning it exactly 1 flat backward.
20. Position the column in the transfer line. Use the septum as a guide to measure the correct length you should insert the column. Be careful not to change the location of the septum on the column.

**Figure 395.** Positioning the Column in the Transfer Line

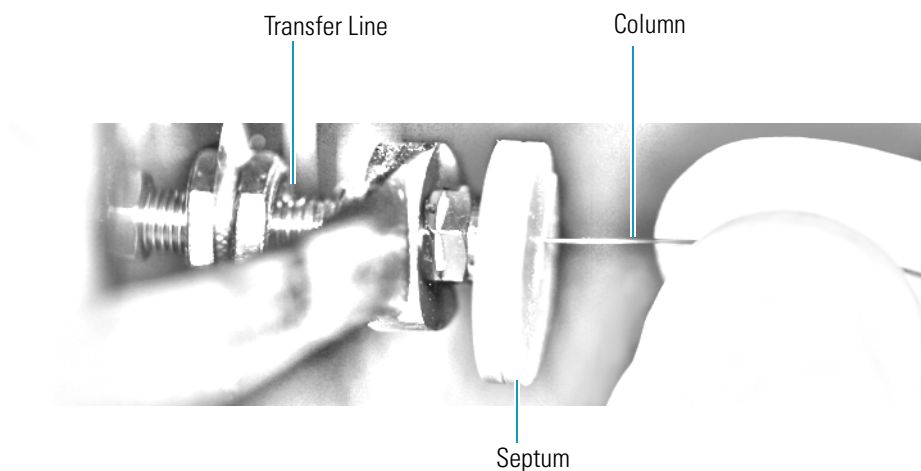

21. Tighten the nut 1 flat forward—back to where it is secure enough in the transfer line that you cannot move it.
22. Tighten the nut 1 additional quarter turn.
23. Remove the cut septum.
24. Close the front door of the GC.

**Note** If you are using a SilTite ferrule, follow the instructions that come with SilTite ferrules. If you are using a graphite Vespel ferrule, they require conditioning to ensure a leak-tight seal. See the *ISQ Spare Parts Guide* for information about ordering these ferrules.

25. Condition the graphite Vespel ferrule:
  - a. Raise the oven temperature to the maximum temperature you will operate the GC column.
  - b. Wait 10 minutes.
  - c. Lower the oven temperature to 40 °C (104 °F) and allow it to cool before continuing.



**WARNING BRUN HAZARD** The oven may be hot. Allow it to cool to room temperature before opening it. The injector will still be hot, so do not touch it.

- d. Retighten the transfer line nut.
26. Restore working conditions.
  - a. Raise the oven temperature to the initial temperature that you will use.
  - b. Turn on vacuum compensation on the GC.
  - c. Twist the vent valve clockwise to close the valve. Be sure not to pin. the o-ring.

- d. If you are using hydrogen as a carrier gas, replace the front panel screw.
- e. Replace all remaining hydrogen safety screws if you are using hydrogen.

27. Power on the mass spectrometer.



**WARNING - FIRE HAZARD:** If you are using hydrogen, do NOT reach over the top of the instrument to power it on. Instead, reach around the side or go to the back of the instrument and flip up the power switch.

28. Once the instrument is pumped down and able to scan, click **Air & Water / Tune** on the Dashboard view air water spectra and look for evidence of leaks with a large  $m/z$  28 signal. If you observe a leak, stop scanning and gently tighten the nut in small increments until no leaks appear when scanning.


## Connecting the Capillaries to the Microfluidics Splitter

You can now attach all the capillaries to the microfluidics splitter.

### ❖ To connect the capillaries to the microfluidics splitter

1. Place the ferrule connected to the original GC column into capillary B location on the microfluidics splitter. See [Figure 396](#). This corresponds to the bottom of the three holes on a correctly installed splitter.

**Figure 396.** Correct Orientation for the Microfluidics Splitter

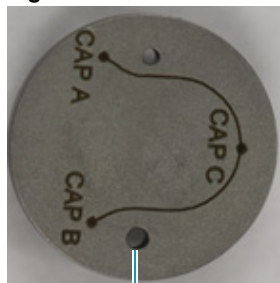


GC Column

**Tip** If a ferrule gets stuck in the microfluidics splitter, use a thumbtack or similar pointed tool for removal. Insert the point of the tool between the ferrule and threads and press so that the ferrule is forced off center. This will dislodge the ferrule.

2. Use the FingerTite tool to secure the nut you previously attached to the column to top of the three holes on the microfluidics splitter.
3. Orient the capillary connected to the transfer line and the one attached to the no-vent module as shown in [Figure 397](#).

**Figure 397.** Correct Orientation of Columns to Microfluidics Splitter

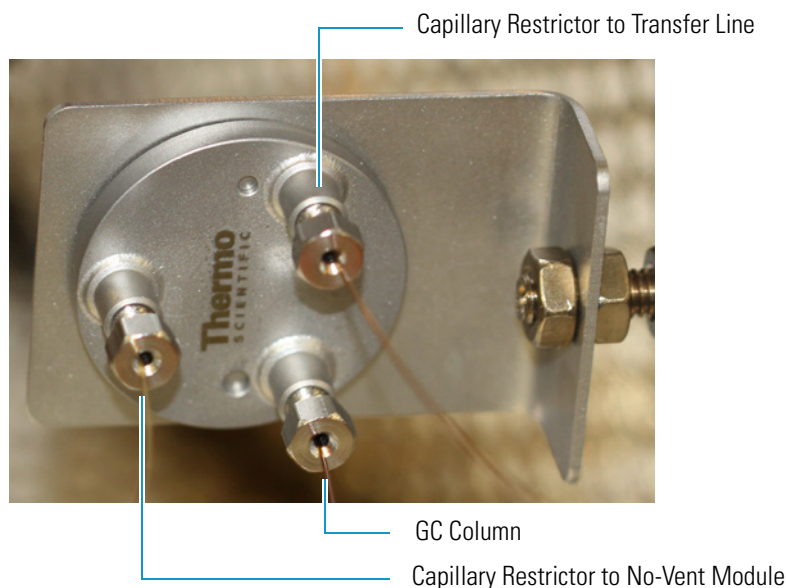

4. Close the GC door.
5. Close the vent valve knob.
6. Power on the mass spectrometer.
7. Ensure the switch is set to **Off** on the NoVent microfluidics module.
8. Let the mass spectrometer pump down for a minimum of one hour and then check the air-water spectrum for gross leaks. Assuming the convection gauge and ion gauge (if present) indicate appropriate pressures, small leaks can be located by spraying with Freon, argon, or another suitable gas near the tubing connections.

## Configuring the Post-Column

After installing the NoVent microfluidics, you must enter the post-column length and ID before running samples. This section contains instructions for configuring your post-column settings on both a TRACE 1310 and a TRACE 1300.

### ❖ To enter post-column settings on a TRACE 1310

1. On the home page of the TRACE 1310 touchscreen, click the **Configuration** icon.

**Figure 398.** Locating the Configuration Icon on the TRACE 1310 Touchscreen



2. The default **GC Configuration (Oven)** screen opens. Click the **Back Column** or **Front Column** button, corresponding to where you installed the No-vent microfluidics.
3. The **GC Configuration (Front or Back Column)** screen opens. Click the down arrow to access the post column settings.
4. Select **Yes** from the **Post Column?** menu. Enter the post column length and ID in the boxes below.

**Figure 399.** Configuring the Post Column Settings




**IMPORTANT** For hydrogen, the **Post Column Length** is 1.2 m. For helium, the **Post Column Length** is 0.6 m. The **Post Column ID** is 0.17 mm for both hydrogen and helium.

5. Click the home icon when you are finished.

❖ **To enter the post-column settings on a TRACE 1300**


1. Open *Xcalibur* by clicking the icon on your desktop
2. The Xcalibur Roadmap opens. Select **TRACE 1300 Series GC** from the list of instruments.
3. Click **Column Setup** in the Status pane.

**Figure 400.** Xcalibur Roadmap

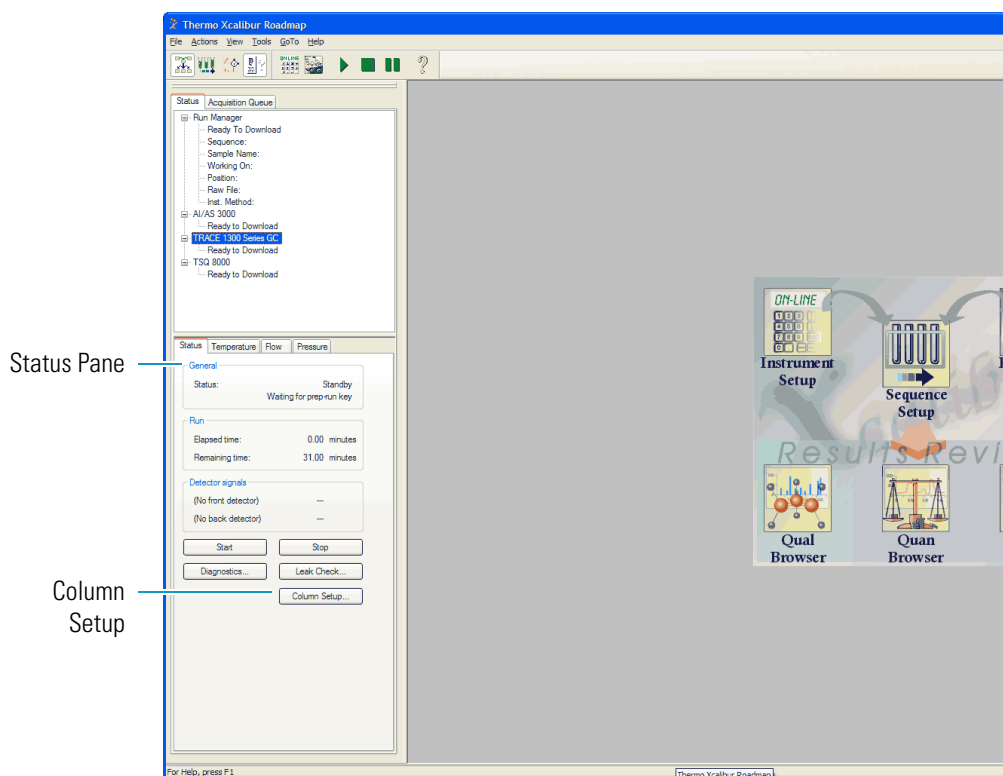

4. The **Column Setup** window opens. Check the **Using Post-Column** checkbox and enter the post-column length and ID.

**Figure 401.** Column Setup Options

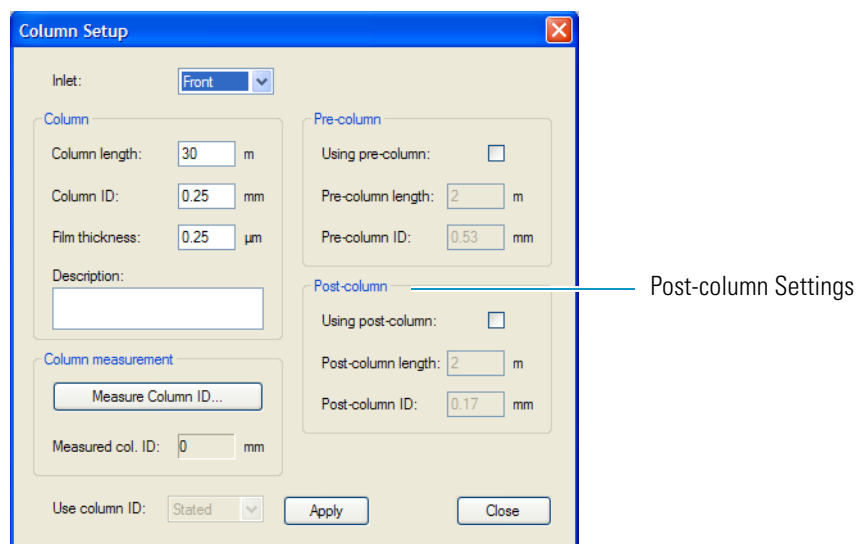


5. Click **Apply**.
6. Click **Close**.

## Using the Module

The GC does not recognize the presence of the NoVent microfluidics, so it is invisible to the system. The module receives only the voltages for supplying the solenoid valve and the On/Off light from the GC.

The normal condition of the module is **Off**. If the switch is in the **Off** position, the Status LED is **off**. See [Figure 402](#).


**Figure 402.** Module in Off Condition

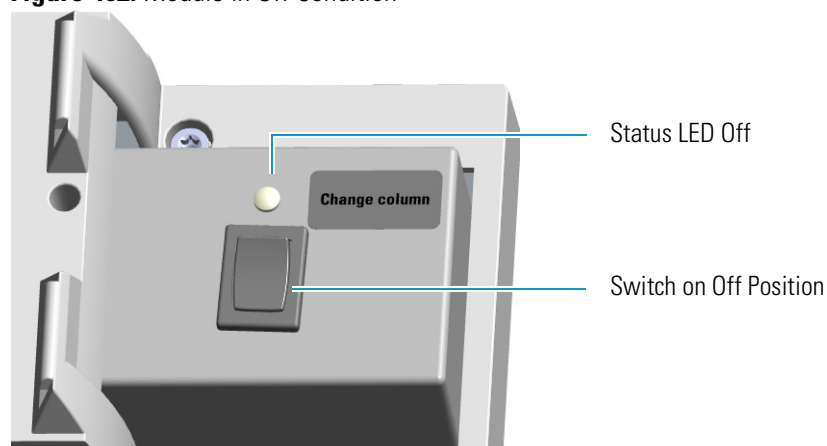

Before replacing the column or performing maintenance on an injector module, flip the switch to the **On** position. Ensure the mass spectrometer is in standby mode and the filament is off. The status LED lights up as solid green indicating that the solenoid valve is activated. The carrier gas flows into the mass spectrometer to avoid breaking vacuum. See [Figure 403](#).



**Figure 403.** Module in On Condition

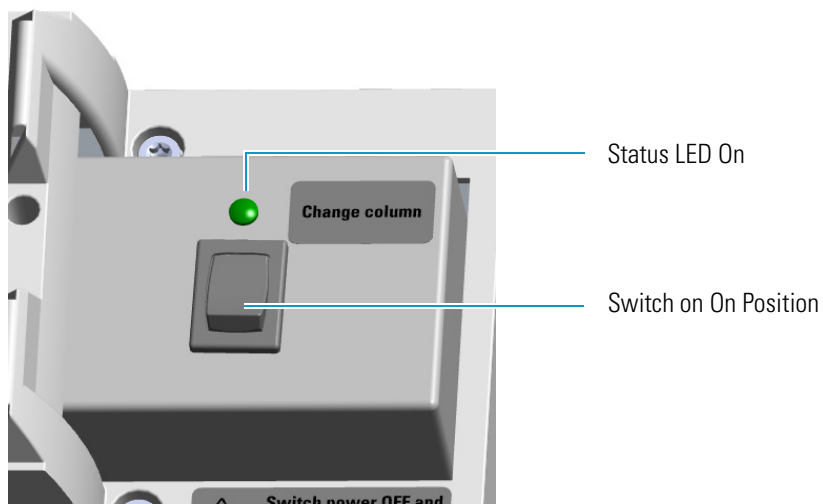


When you have completed your tasks, flip the switch to the **Off** position. The status LED turns **Off** indicating the solenoid valve is deactivated.

## Installing a FTIR Make-up Module

This section provides instructions for installing the FTIR Make-up module on your TRACE 1300/1310 for the connection with the Thermo Scientific Nicolet iS50 FT-IR Spectrometer. According to the configuration of your TRACE 1300/TRACE 1310, a dummy module is present into the free site where the detector module is not installed. See [Figure 404](#).

**Figure 404.** Add a Front/Back Detector Module

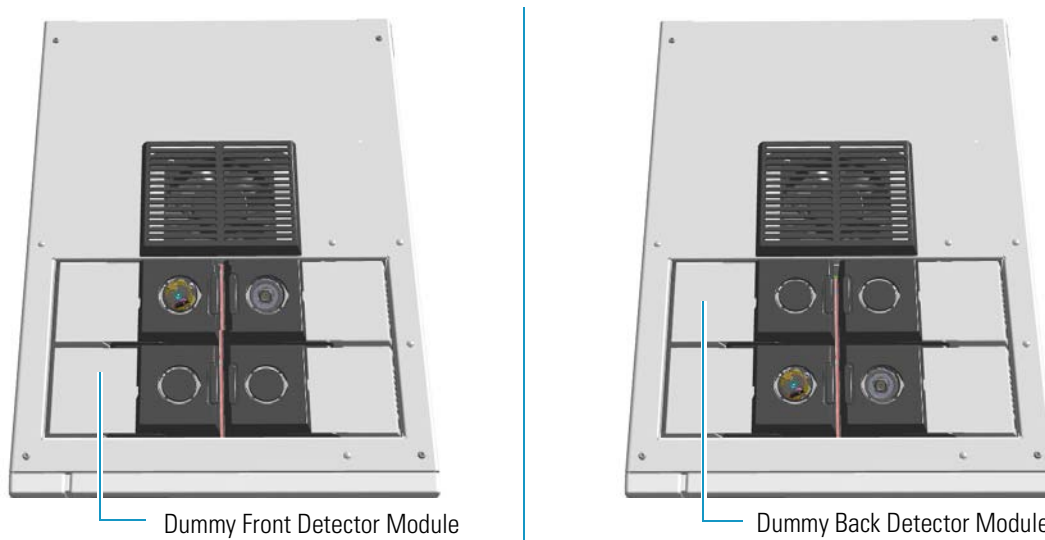

**ATTENTION** Where a dummy module is installed, the gas connection is blocked by a plug.

The FTIR Make-up module must be installed into the detectors housing, and supplied through the line of the front/back carrier gas.

**IMPORTANT** If the front/back carrier gas is not already used by an injector module, you must install the front/back carrier gas inlet, and connect the carrier gas line as described in the section [“Making the Gas Supply Plumbing Connections”](#) on [page 7](#).

The module is equipped with an ON/OFF needle valve for the regulation of the gas flow.

- Open the needle valve manually turning the knob counter-clockwise.
- Close the needle valve manually turning the knob clockwise.

Perform the regulation of the gas flow using a little flat head screwdriver.

- Increase the gas flow turning the flow regulation screw counter-clockwise.
- Decrease the gas flow turning the flow regulation screw clockwise.


The standard operating value is in the range 0.2- 0.3 mL/min. Refer to the *Nicolet iS50 FT-IR Spectrometer Manual* for details.



**ATTENTION** The flow of the carrier gas depends on the inlet pressure. Set the carrier gas flow according to the working pressure.

The components of the FTIR module are shown in [Figure 405](#).

**Figure 405.** FTIR Make-up Module

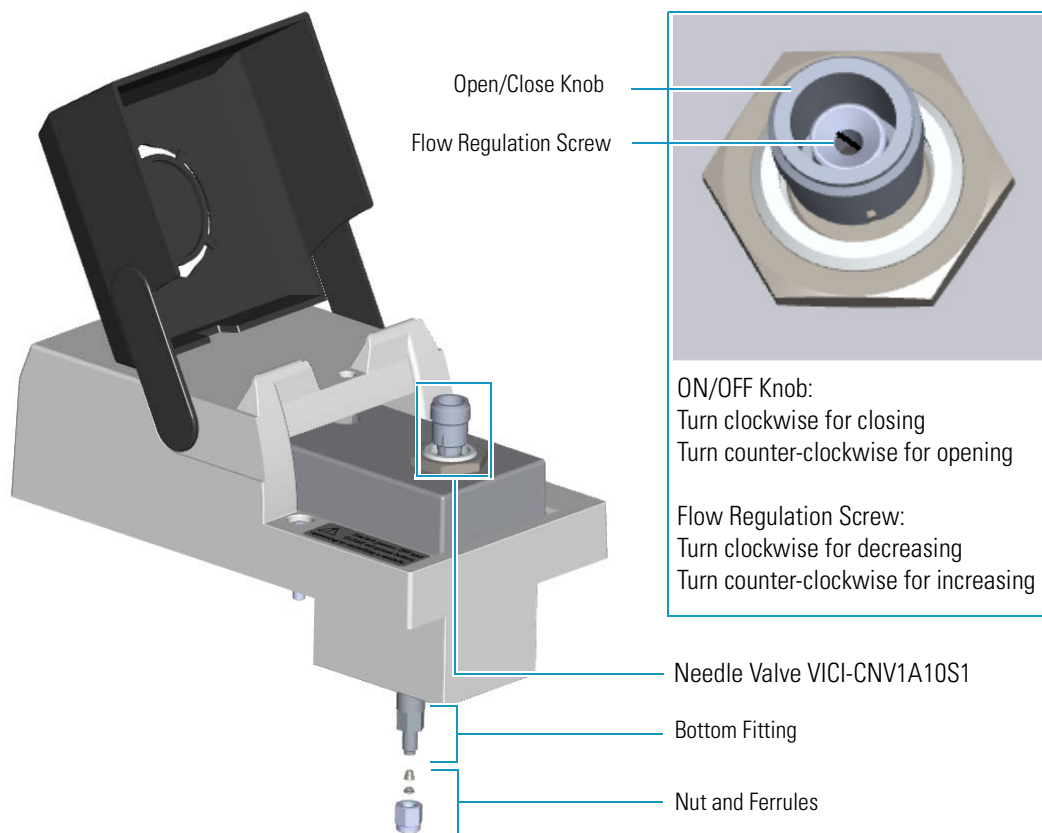

❖ **To add a FTIR Make-up module**


1. Put the GC in standby condition.
2. Cool the oven, injector, and detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Close the gas supplies.
4. Power off the GC.

- a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
6. Remove the dummy module from the position where the FTIR module will be installed.
- a. Open the module flap cover.
  - b. Using a T20 Torxhead screwdriver, unscrew and remove the two captive fixing screws.
  - c. Keeping the dummy module flap cover open, lift up the module from its seat in the detector housing. Place the dummy module on a clean surface.
  - d. Remove the gas block plug from the gas connections by unscrewing its fixing screw using a T20 Torxhead screwdriver. See [Figure 406](#).

**Figure 406.** Detector Gas Block Plug

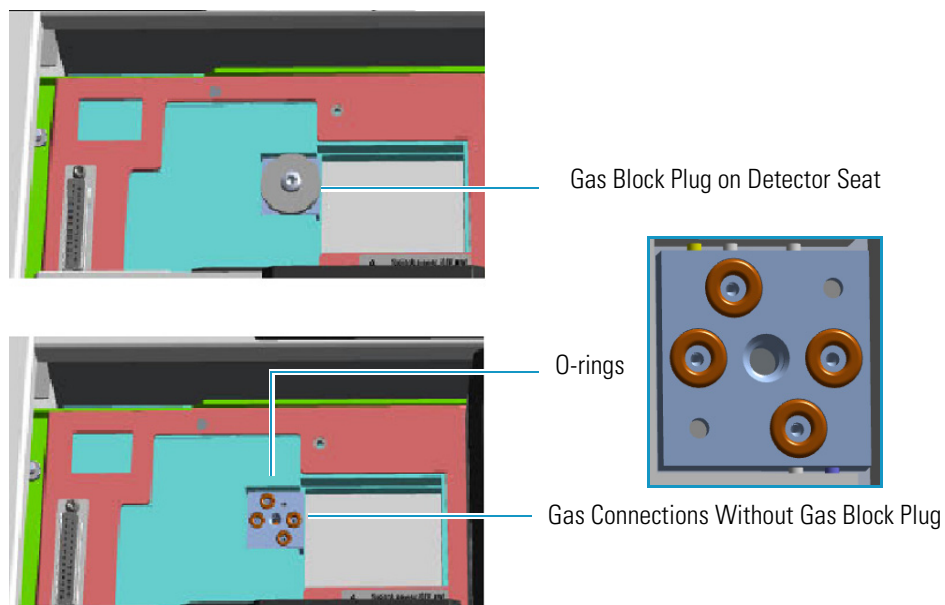

**WARNING** Make sure all the four o-rings are placed into their seats on the gas connection. See [Figure 406](#).

**Do not install the module if the o-rings are missing.**

7. Open the front door of the GC.
8. Plug the FTIR Make-up module into the main frame.
  - a. Open the module flap cover.

- b. Keeping the module flap cover open, place the module in its seat.
- c. Use a T20 Torxhead screwdriver to tighten the three captive fixing screws without overtightening.




**ATTENTION** To maintain the correct alignment the screws must be tightened in turn. Tighten each screw only a small amount before moving to the next screw. Repeat until all are secure.

- d. Close the module flap cover.
9. Connect the 1/16-in tubing provided.
- a. Connect an end of the 1/16-in tubing to the bottom fitting of the module using the Swagelock® nut and ferrules provided.
  - b. Connect the other end of the tubing to the transfer line of the FT-IR spectrometer. Refer to the *Nicolet iS50 FT-IR Spectrometer Manual* for details.
10. Open the gas supplies.
11. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
12. Power on the GC.
- a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
13. Close the front door of the GC.

# Installing the Hot Injection Adapter Kit on the SSL/SSLBKF Injector


This section provides the instruction for installing the Hot Injection Adapter kit (PN 19050733) on the SSL/SSLBKF injector for the use with an external device for gas sampling. See [Figure 330](#).

**Figure 407.** Hot Injection Adapter

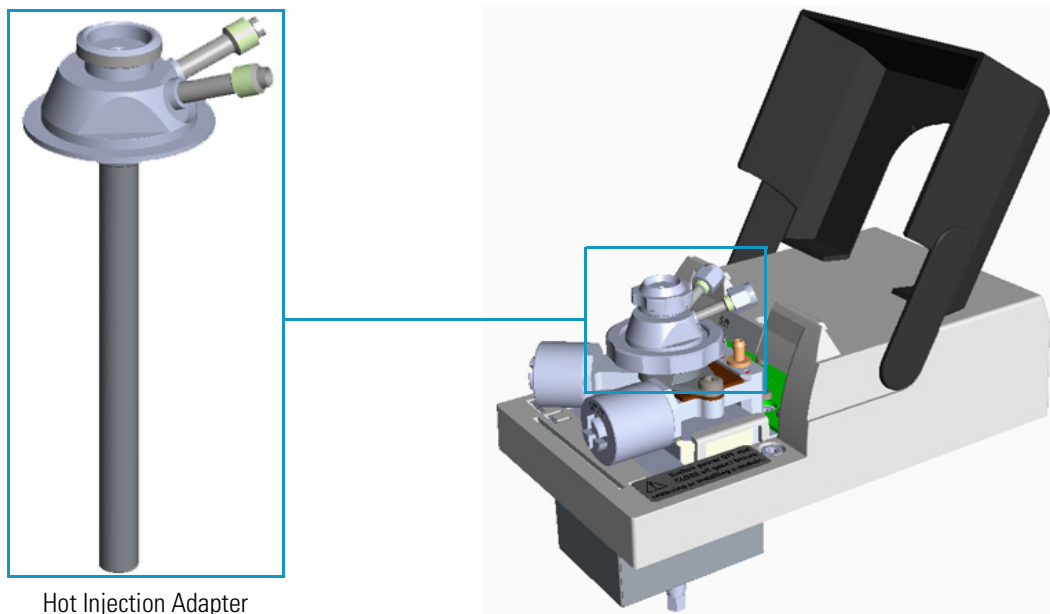

The Hot Injection Adapter consists of the following parts: head, short PTV liner, liner seal, and stem. See [Figure 408](#).


**Figure 408.** Hot Injector Adapter Assembly

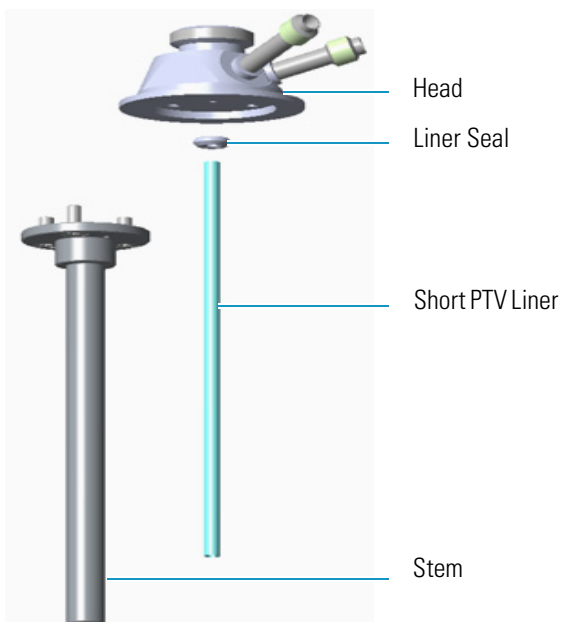

❖ **To install the hot injection adapter on the injector**

1. Assemble the hot injection adapter if not already done. See [Figure 409](#).

**Figure 409.** Adapter Assembly

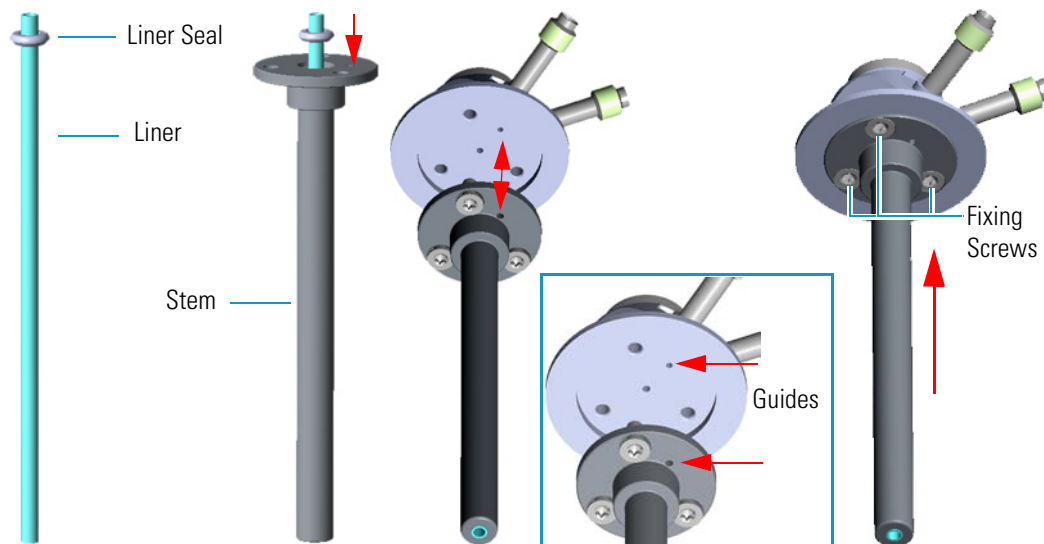

- a. Place the liner seal on the liner.
- b. Insert the liner into the stem.
- c. Couple the head and the stem paying attention to properly align the corresponding holes doing match the polarizing guides on the head and on the flange of the stem.
- d. Fix the stem to the head screwing the three fixing screws using a T8 Torxhead screwdriver.

**Note** The stem conducts heat from the body of the SSL/SSLBKF injector to the head of the adapter allowing a hot injection.

2. Put the GC in standby condition.
3. Cool the oven, injector, and detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Close the gas supplies.
5. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
6. Remove the column end from the injector.

## 8 Installing Optional Kits

Installing the Hot Injection Adapter Kit on the SSL/SSLBKF Injector


- a. Open the front door of the GC.
  - b. Loosen the retaining nut from the injector fitting on the upper interior wall of the GC oven.
  - c. Remove the analytical column with its nut and ferrule from the bottom of the injector.
7. Put the autosampler away if present.
  8. Remove the SSL/SSLBKF injector module from its seat.
    - a. Open the module flap cover.
    - b. Using a T20 Torxhead screwdriver, unscrew the three captive fixing screws.
    - c. Throw upward the module from its seat of the injector housing. Place the SSL/SSLBKF module on a clean surface.

**WARNING** Make sure the o-ring remains into its seat on the gas connection.  
**Do not install the module if the o-ring is missing.**



9. Replace the module flap cover with the one provided. See [Figure 331](#).

**Figure 410.** Module Flap Cover for Hot Injection

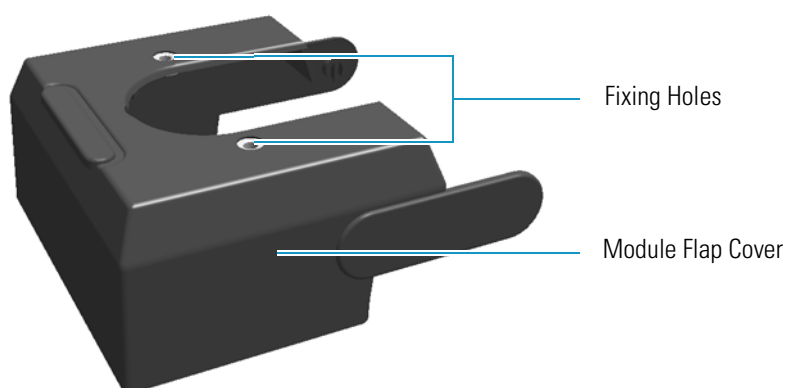

**Note** The fixing holes present on the top of the module flap cover can be used for making more stable a sampling device that requires to be seat directly on the injector.

10. Reinstall the injector module into the main frame.
  - a. Open the module flap cover.




- b. Place the injector module in its seat. Be sure to insert the 25-pin male connector, on the bottom of the module, into the 25-pin female connector on the injector seat of the injector housing.
- c. Use a T20 Torxhead screwdriver to tighten the three captive fixing screws without overtighten.



**CAUTION** To maintain the correct alignment the screws must be tightened in turn. Tighten each screw only a small amount before moving to the next screw. Repeat until all are secure.

11. Remove the top parts of the injector. See [Figure 332](#).

**Figure 411.** Injector Top Parts Removal

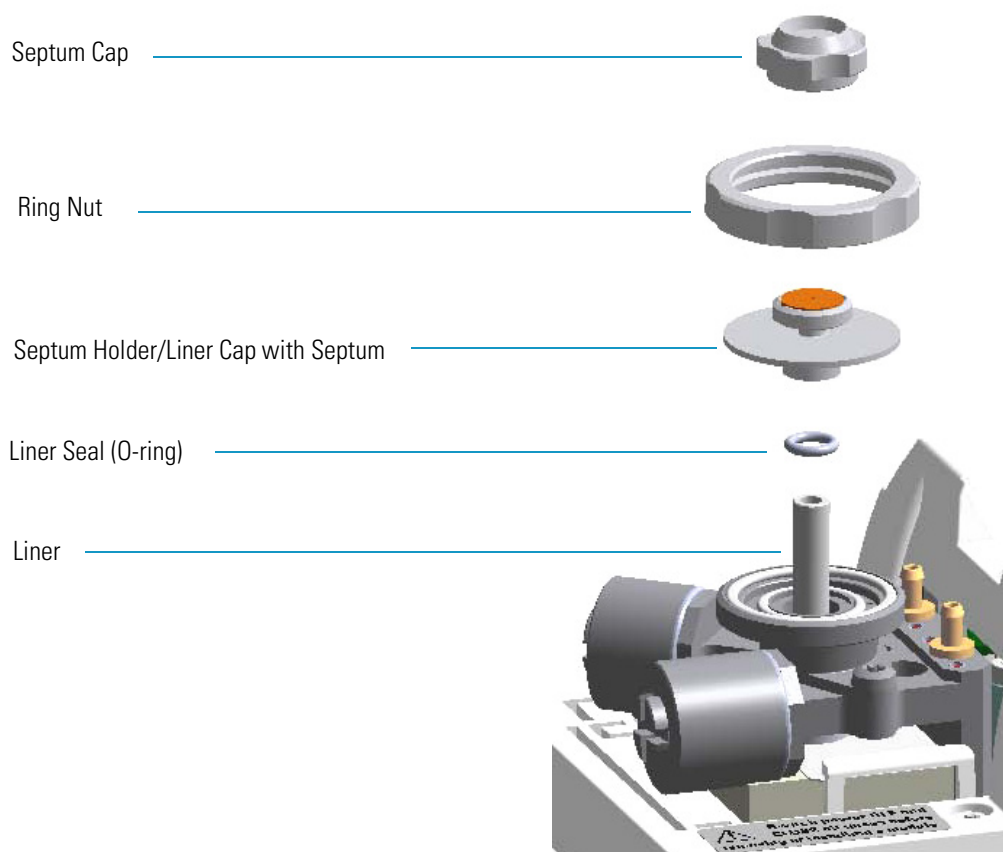


- a. Unscrew the septum cap of the injector.
  - b. Unscrew the ring nut.
  - c. Remove the septum holder/liner cap with septum from the injector body head.
12. Remove the liner and its seal.
    - a. Use tweezers to remove the current liner with the liner seal (o-ring) from the injector.
  13. Remove the injector body. See [Figure 412](#).

## 8 Installing Optional Kits

Installing the Hot Injection Adapter Kit on the SSL/SSLBKF Injector

**Figure 412.** SSL Injector Body Removal

SSL/SSLBKF Injector Body Fixing Screws

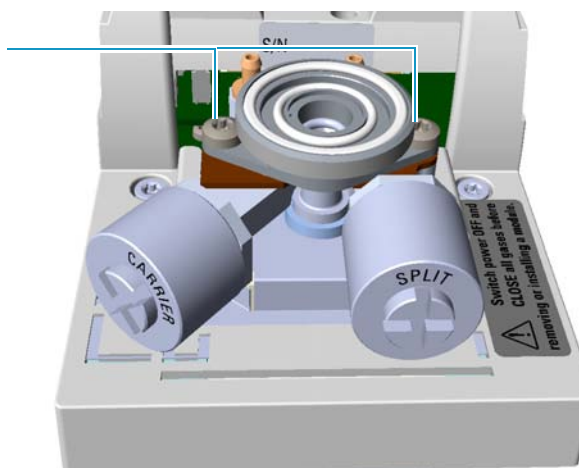

External Body Head O-ring


Internal Body Head O-ring

SSL Injector Body

Carrier Line O-ring

Purge Line O-ring


Split Line O-ring

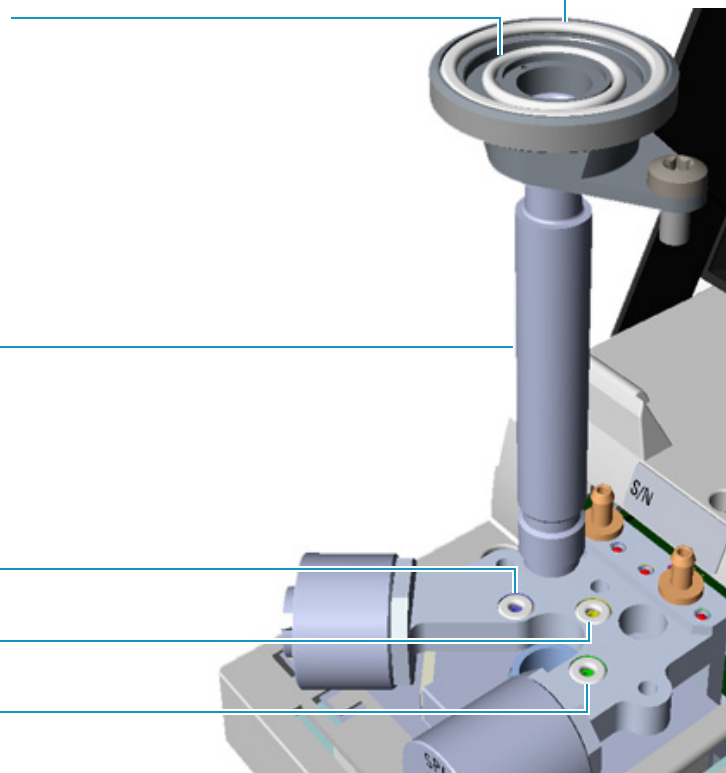

- a. Using a T20 Torxhead screwdriver, undo the two injector body fixing screws, and extract the injector body from its housing.

**Note** Do not remove the carrier, split, and purge lines o-rings.


14. Install the spacer plate. See [Figure 413](#).

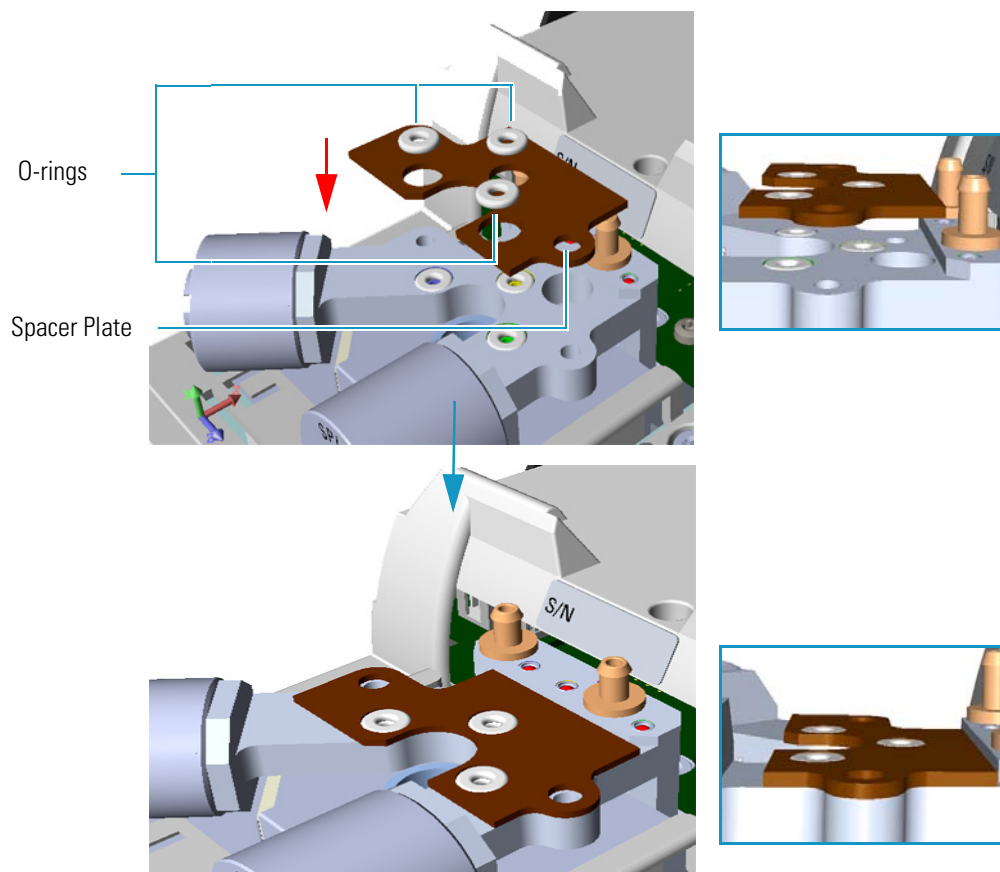**Figure 413.** Spacer Plate

- a. Place the spacer plate and its three o-ring on the injector base.

**Note** The three o-ring of the spacer plate do not replace the carrier, purge, and split lines o-rings but are put on them, then six o-rings are present on the base.

15. Reinstall the injector body. See [Figure 414](#).

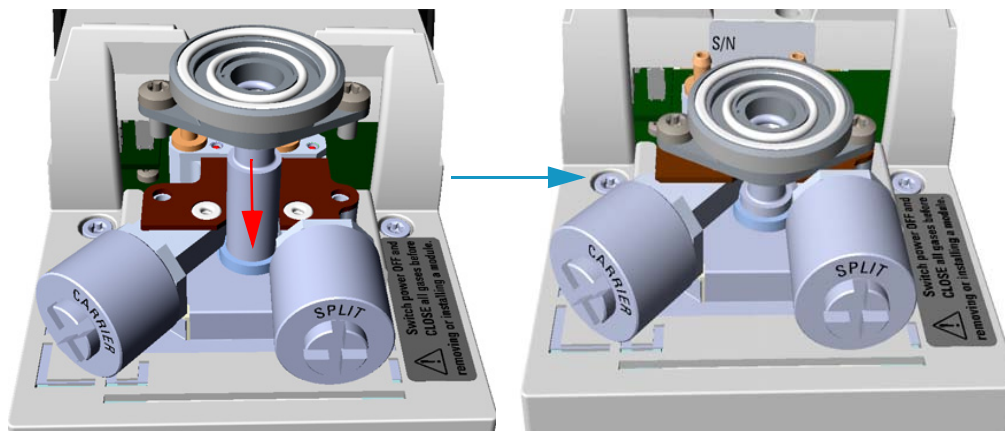**Figure 414.** Injector Body Reinstallation


- a. Reinstall and fix the injector body into its housing by screwing the two fixing screws.

## 8 Installing Optional Kits

Installing the Hot Injection Adapter Kit on the SSL/SSLBKF Injector


- b. If required, replace both the internal and external body head o-rings using tweezers.
- 16. Install the hot injection adapter on the top of the injector. See [Figure 333](#).

**Figure 415.** Hot Injection Adapter (1)

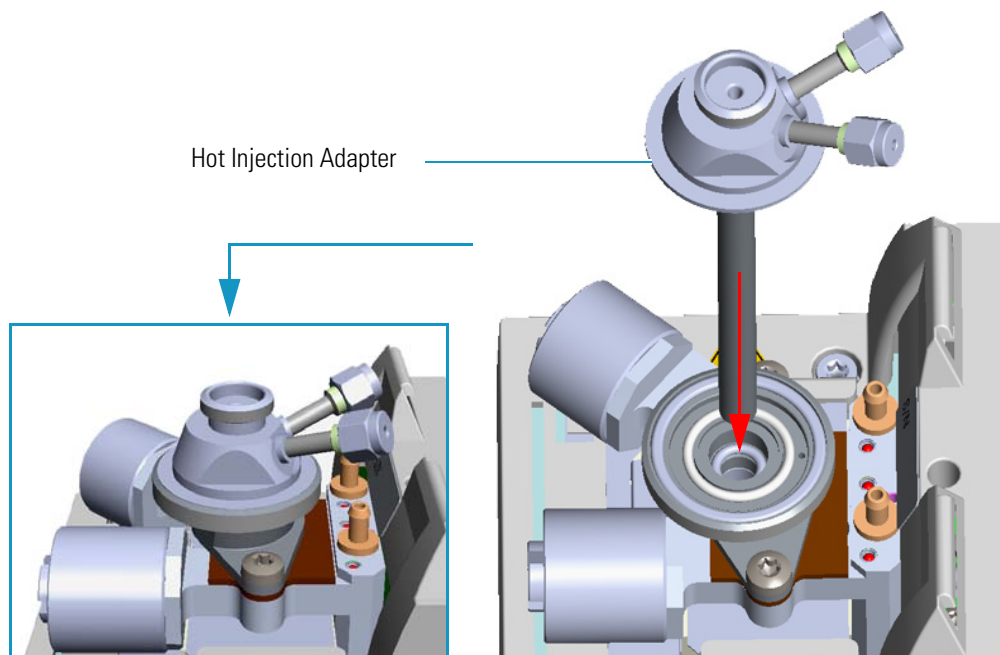


- a. Avoid touching the septum with your fingers. Insert a new septum into the septum holder cavity of the hot injection adapter using tweezers. See [Figure 334](#).


**Figure 416.** Hot Injection Adapter (2)

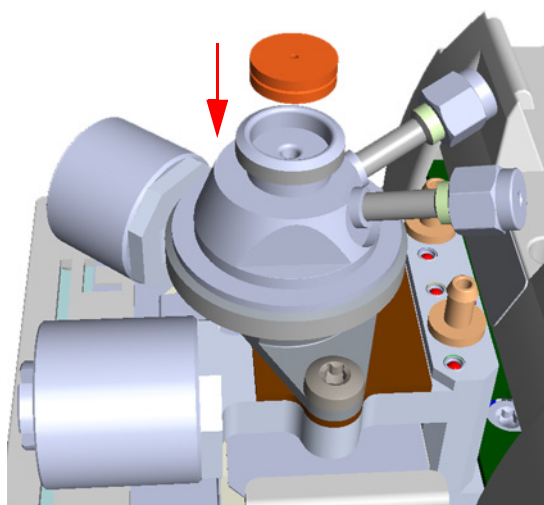

- b. Guide the ring nut on the hot injection adapter with the septum, then fix it screwing the ring nut. See [Figure 336](#).

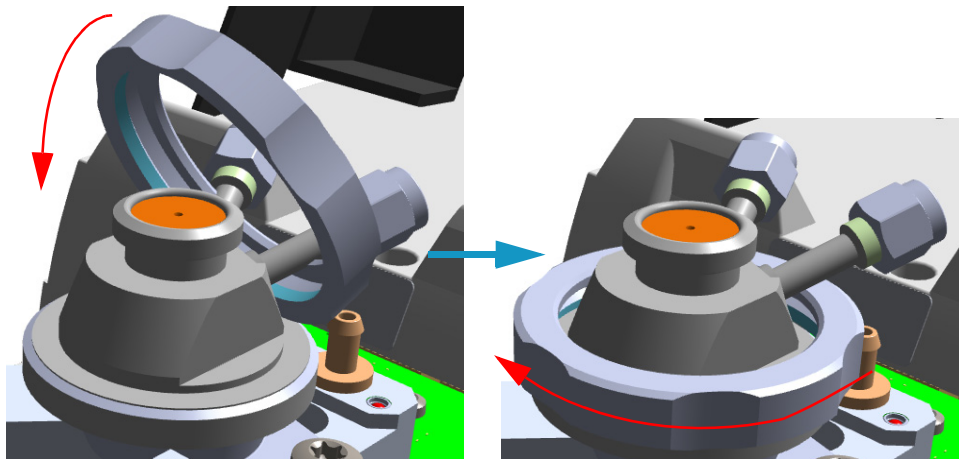**Figure 417.** Hot Injection Adapter (3)

- c. Screw and tighten the septum cap to finger-tight. See [Figure 418](#).

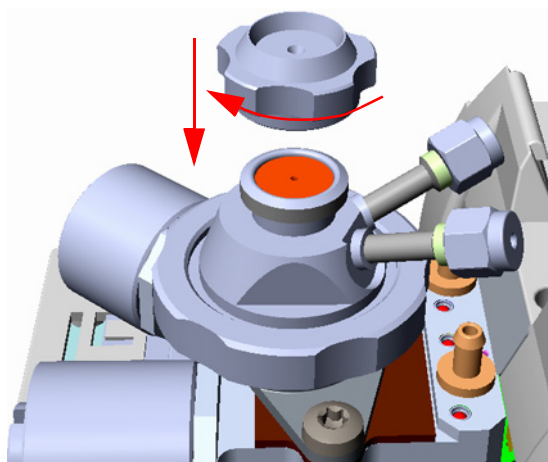**Figure 418.** Hot Injection Adapter (4)


**CAUTION** Do not overtighten the septum cap because you might damage the septum and affect instrument performance.

17. Close the module flap cover.
18. Connect the hot injection adapter to the device for the sampling of gas. See the example in [Figure 419](#).

## 8 Installing Optional Kits

Installing the Hot Injection Adapter Kit on the SSL/SSLBKF Injector

**Figure 419.** Hot Injection Adapter (5)

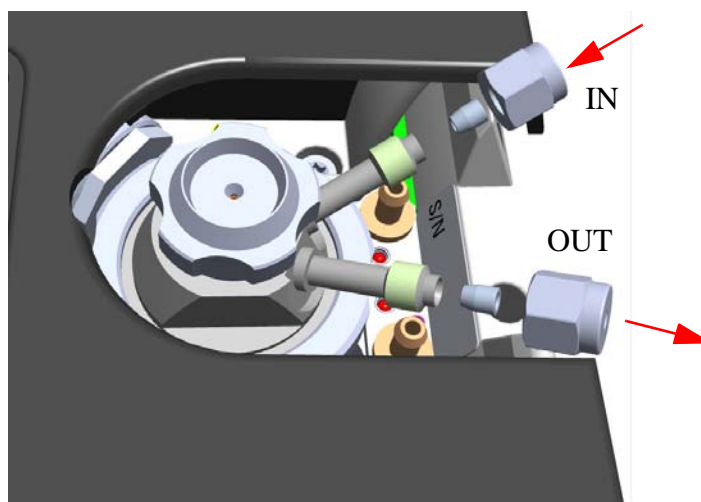


**Note** If only a port of the adapter is connected to the device for the sampling of gas, close the other one with the blind nut provided.

19. Reconnect the column end to the injector and verify the connection point.
20. Open the gas supplies.
21. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip out the power switch (breaker) to the position I.
22. Pressurize the module with the carrier gas.
23. Check for leaks.
  - a. Use a handheld electronic leak detector (Thermo Scientific GLD Pro leak detector or equivalent) to check the two fitting for leaks.
  - b. If you detect a leak, tighten the connection and retest it.
  - c. Repeat this process until all connections are leak free.
24. Close the front door of the GC.
25. If present, update the autosampler for the new injection position.

## Performing the Dual FPD Detector Configuration

This section provides the instructions for performing the Dual FPD configuration by connecting a second photomultiplier tube with a different interferential filter on the same detector control module. See [Figure 420](#).

**Figure 420.** Dual FPD Configuration

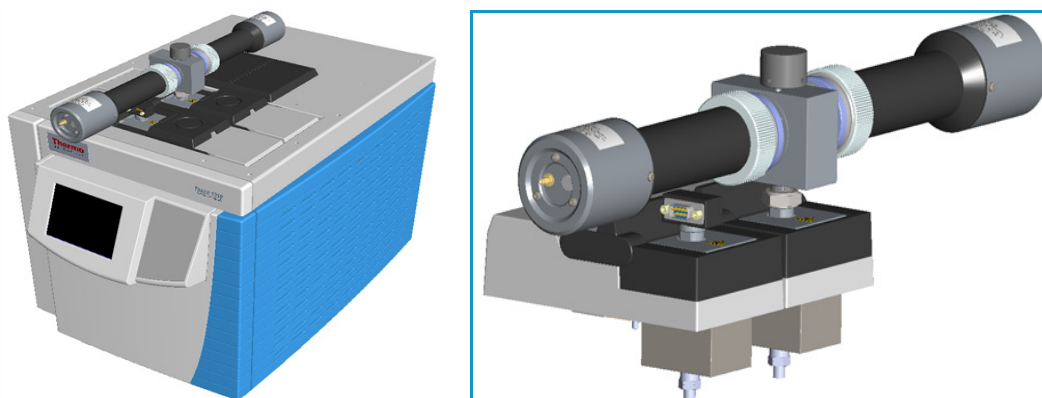

To perform the Dual FPD configuration you need the following devices:

- Two Flame Photometric Detector modules
- A Dual FPD upgrade kit PN 19050783



**IMPORTANT** The Dual FPD must be installed onto the FPD control module placed into the **back** site of the detectors housing.

The FPD control module placed into the **front** site of the detector housing must be equipped with the **PT100 simulator**, provided with the kit, instead of the FPD detector, and will supply the second photomultiplier of the configuration.

The operation consists of the following steps:

- Remotion of the **front** and **back** FPD Detectors from the relevant detector control module.
- Disassembling the photomultiplier tube from a FPD detector.
- Assembling the second photomultiplier tube on the other FPD detector to obtain the Dual FPD configuration.
- Remotion of the **jet** from the **front** detector control module.
- Installation of the **PT100 simulator**, provided with the kit, on the detector base of the front detector control module.



**Note** For writing convention, we consider the disassembling of the photomultiplier tube from the FPD detector installed on the **front** detector control module, and the assembling of the Dual FPD configuration on the FPD detector installed on the **back** detector control module. We name them **secondary FPD detector** and **primary FPD detector** respectively.

#### ❖ To assembly the Dual FPD configuration

##### Materials needed

Dual FPD Upgrade Kit PN 19050783

Gloves

FPD fixing tools

Appropriate interferential filter (394 nm for Sulphur; 526 nm for Phosphorous; 610 nm for Tin

Methylene chloride or GC-grade methanol

Screwdriver

T6 Torxhead key

T10 Torxhead screwdriver

1-mm Allen wrench

5 mm wrench

5/16-in tube wrench

#### Initial operations

1. Put the GC in standby condition.
2. Switch off the flame. Hydrogen and air are automatically closed.
3. Cool the detector to room temperature.


**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If other external modules are present, unplug the power cable from the AC Input connector of each module, and from the wall outlet.
6. Remove the FPD detectors from the base of the respective detector control module.




- a. Disconnect the signal, excitation voltage, and ignition/heating cables from the detector. See [Figure 421](#).

**Figure 421.** FPD Cable Disconnection

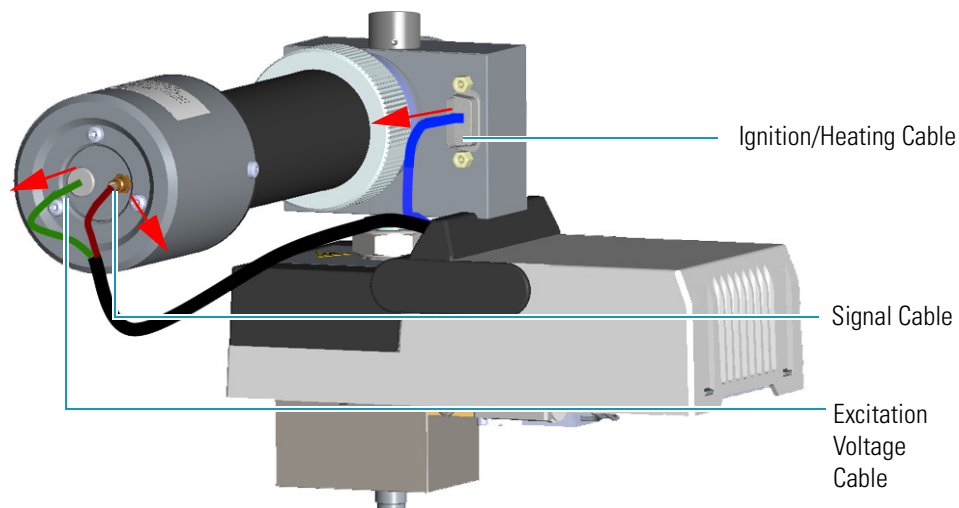

- b. Using the tool provided with the system, loosen the fixing nut on the base of the detector and remove it. See [Figure 422](#).

**Figure 422.** FPD Detector Removal

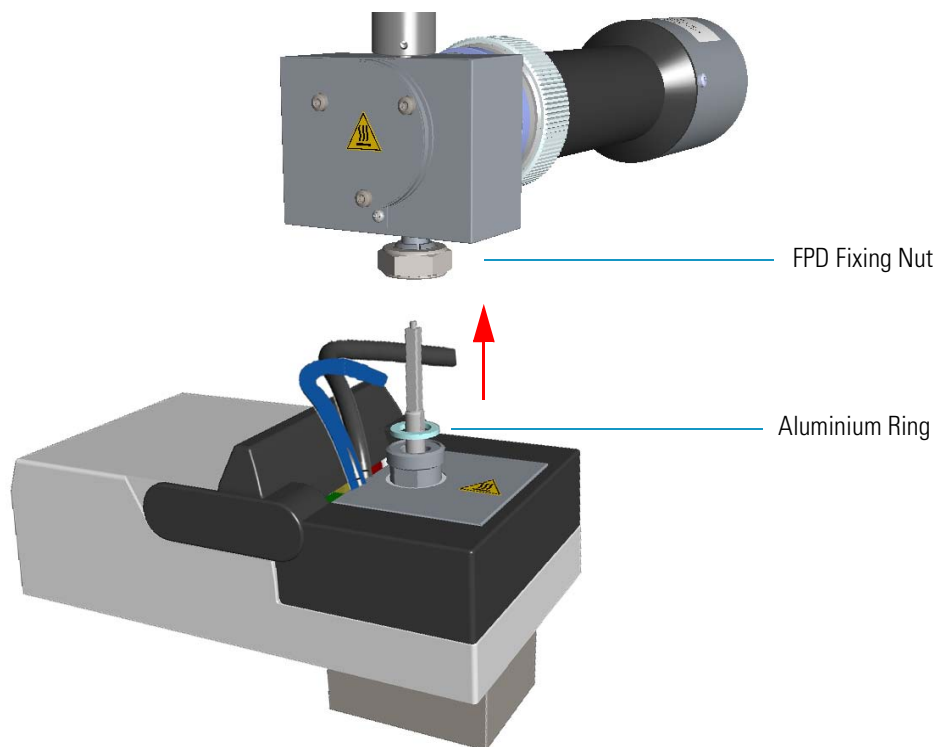


**Note** Do not lose the aluminium ring inserted between the detector head and the base body.

- c. Place both the FPD detectors on a clean and free workbench.

#### Disassembling the secondary FPD detector

7. Remove the photomultiplier tube assembly and the filter from the cell of the secondary FPD detector. See [Figure 423](#).

**Figure 423.** Photomultiplier Assembly Removal

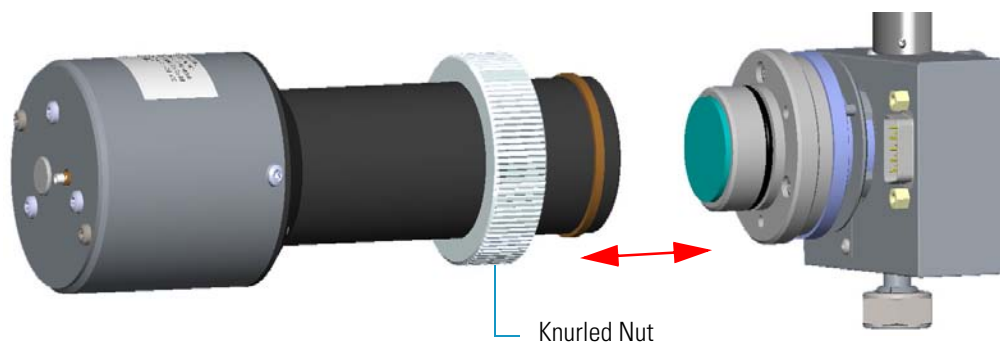


- a. Loosen the knurled nut that fixes the photomultiplier assembly and remove it from the detector body.



**CAUTION** The photomultiplier tube could damage if exposed to ambient light with the excitation voltage On. Make sure the power supply has been switched off before disconnecting the tube from the detector body.


- b. Remove the interferential filter from its housing, handling it very gently. Keep it using a clean paper towel. See [Figure 424](#).

**Figure 424.** Filter Removal

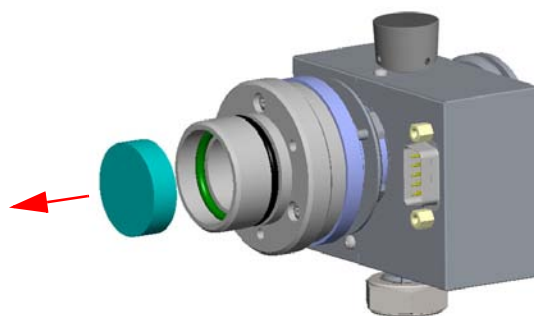


**CAUTION** Filters are fragile. Pay attention not to let the filter fall down and damage.

**CAUTION** If the filter installed is equipped with the spacer (see the kits PN 19050785 and PN 19050786), first remove the spacer, and then the filter.

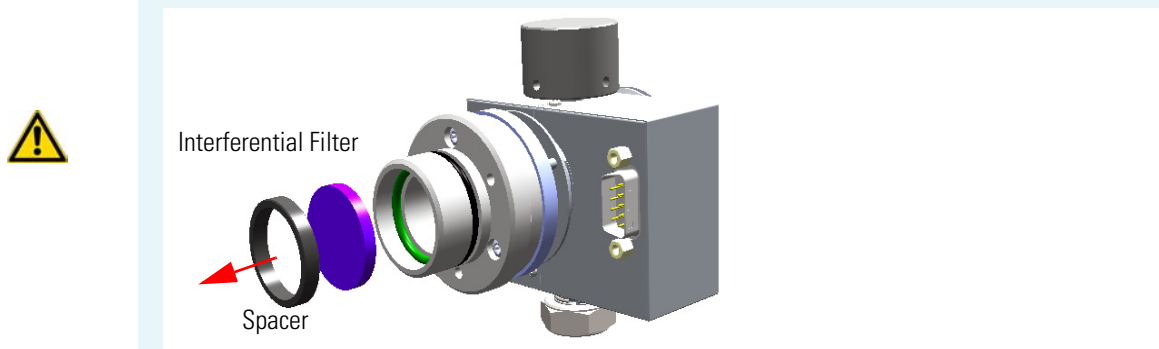


8. Remove the filter a support and the flange assembly from the cell of the secondary FPD detector.
  - a. Loosen the three smaller Allen screws that fix the flange to the detector body and remove it with the spacer that should remain inserted in the flange See [Figure 425](#).

**Figure 425.** Filter Support and Flange Assembly Removal

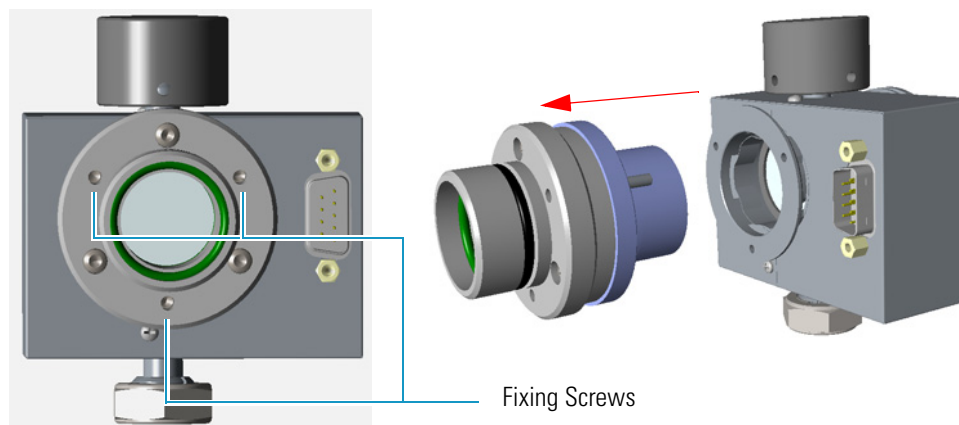

- b. Now you could access the flame-side heat shield, and the relevant graphite seal. Loosen the Allen screws that fix the mirror plug and remove it. See [Figure 426](#).

**Figure 426.** Mirror Plug Removal

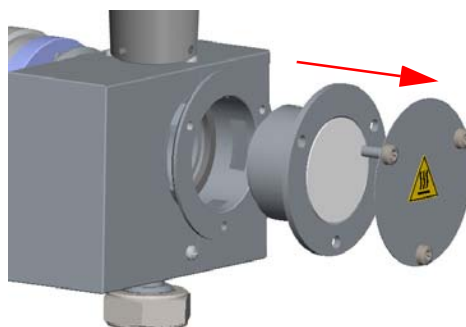


- c. Insert the handle of a screwdriver or other un-sharpened tool in the combustion chamber and push the flame-side heat shield and its graphite seal out from its housing. Act gently to avoid breaking of the heat shield. See [Figure 427](#).



**CAUTION** While pushing out the heat shield, pay attention not to damage the ignition coil.

**Figure 427.** Filter-side Heat Shield Removal

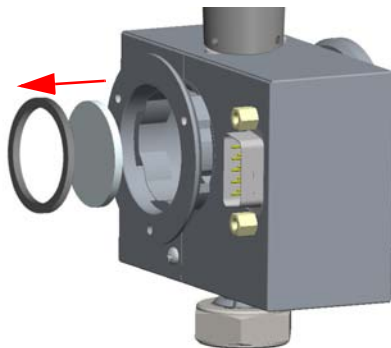


- d. Remove with care traces of graphite due to the breaking of the seal.
- e. Reinstall the mirror plug.

#### Assembling the primary FPD detector in Dual FPD configuration

Perform the Dual FPD configuration on the primary FPD detector using the components removed from the secondary FPD detector.

9. Remove the mirror plug from the cell of the primary FPD detector.
  - a. Loosen the Allen screws that fix the mirror plug and remove it. See [Figure 428](#).


**Figure 428.** Mirror Plug Removal

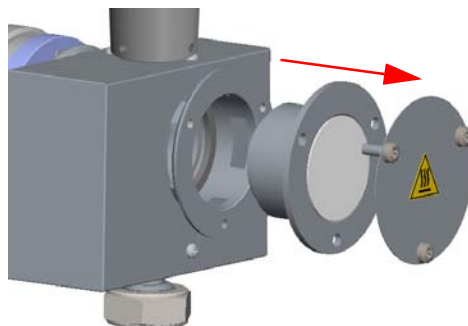

10. Insert the filter support and the flange assembly into the primary FPD detector.
  - a. Insert the heat shield into its housing inside the detector body.
  - b. Insert the spacer into the flange, letting it stand out for about 5 mm.

- c. Insert the graphite seal on the spacer, pushing it slowly until it touches the flange.
- d. Fix the assembly tightening the three Allen screws. See [Figure 429](#).


**Figure 429.** Filter Support and Flange Assembly Installation

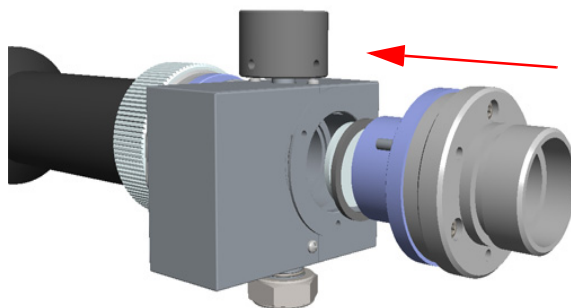

11. Place the appropriate filter and the photomultiplier tube assembly.


- a. Reinsert the filter into the support. The mirror face must be oriented towards the flame. See [Figure 430](#).

**Figure 430.** Filter Remounting

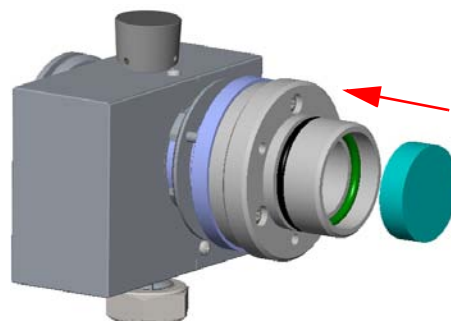

**CAUTION** Avoid touching the filter with your fingers. If you see fingertips on the filter, clean it using a clean paper towel and, if needed, GC-grade methanol before remounting.

**CAUTION** If the filter installed is equipped with the spacer (see the kits PN 19050785 and PN 19050786), first insert the filter, and then the spacer.

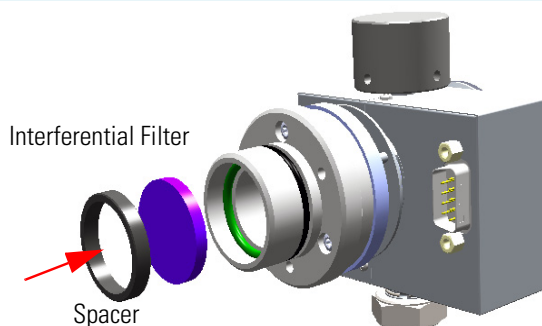


## 8 Installing Optional Kits

### Performing the Dual FPD Detector Configuration


- b. Mount the photomultiplier tube assembly then fix it tightening the knurled nut. See [Figure 431](#).

**Figure 431.** Photomultiplier Tube Assembly Mounting

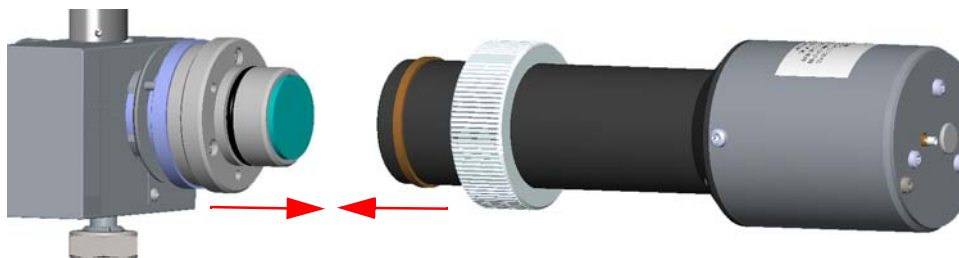

The result of the operation is shown in [Figure 432](#).


**Figure 432.** Result of the Dual FPD Configuration

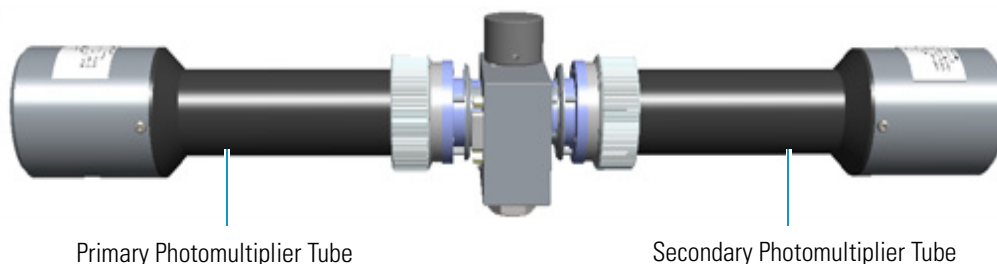

### Installing of the Dual FPD and electrical connections

12. Remove the jet of the front detector control module.
  - a. Using the 5 mm wrench provided with the GC, loosen the jet and remove it from the detector base body. See [Figure 433](#).

**Figure 433.** FPD Jet Removal

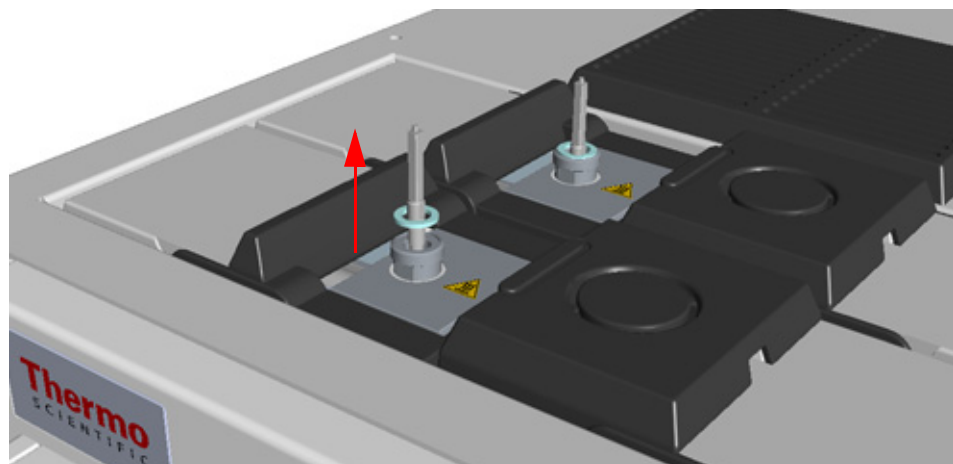




**IMPORTANT** Save the jet and the aluminium ring in a safe place because will be reused when you restore the original configuration of the FPD modules.


13. Install the PT 100 simulator, provided with kit PN 19050783, on the detector base of the front detector control module. See [Figure 434](#).

**Figure 434.** PT100 Simulator Installation

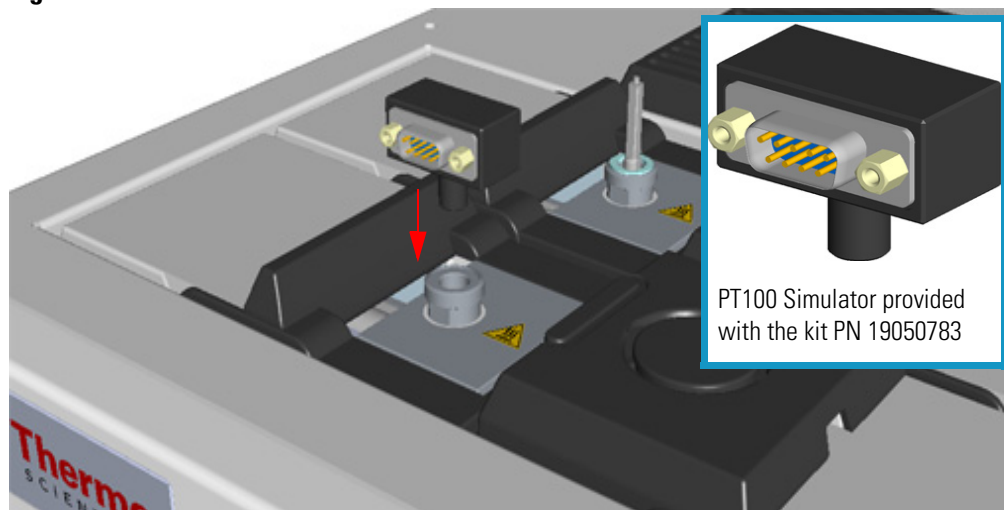

14. Install the Dual FPD on the detector base of the detector control module placed into the back site of the detector housing. See [Figure 435](#).

**Figure 435.** Dual FPD Installation

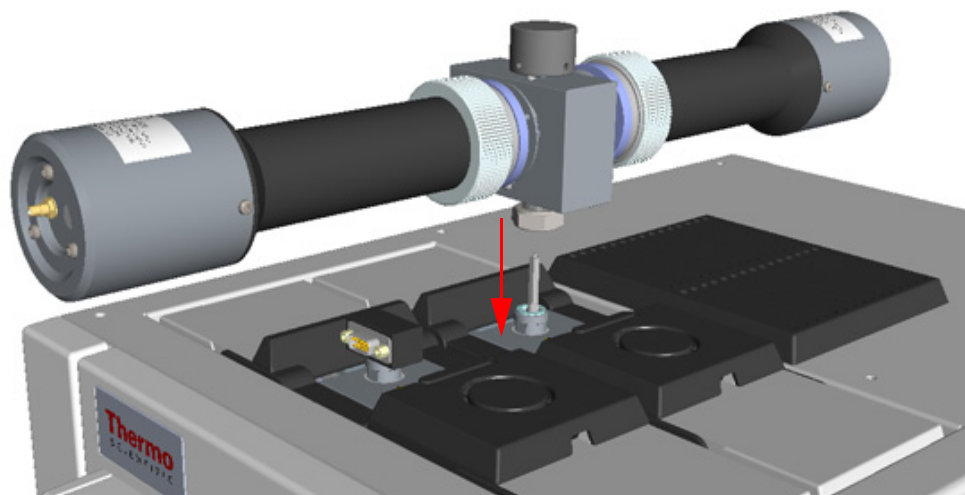


- a. Place the detector on its base body, ensuring that the aluminium ring has been inserted in the correct position, then tighten the fixing nut.
15. Connect the cables to the FPD detectors.

## 8 Installing Optional Kits

### Performing the Dual FPD Detector Configuration

- a. Connect the signal, excitation voltage and ignition/heating cables coming from the control module installed into the back site to the primary FPD detector. See [Figure 436](#).


**Figure 436.** Primary FPD Detector Cables Connection

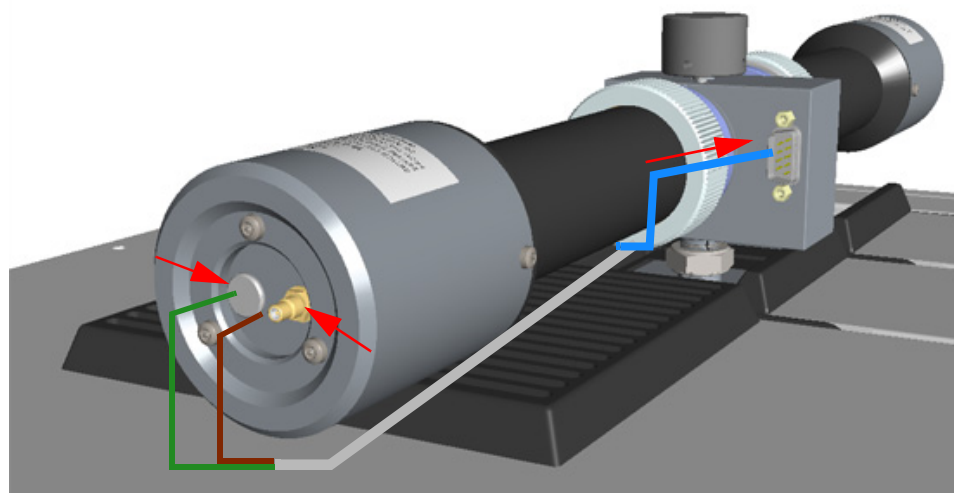

Excitation Voltage Cable   ■   Signal Cable   ■   Ignition Heating Cable   ■

- b. Connect the signal, excitation voltage and ignition/heating cables coming from the control module installed into the front site to the secondary FPD detector. See [Figure 436](#).

**Figure 437.** Secondary FPD Detector Cables Connection

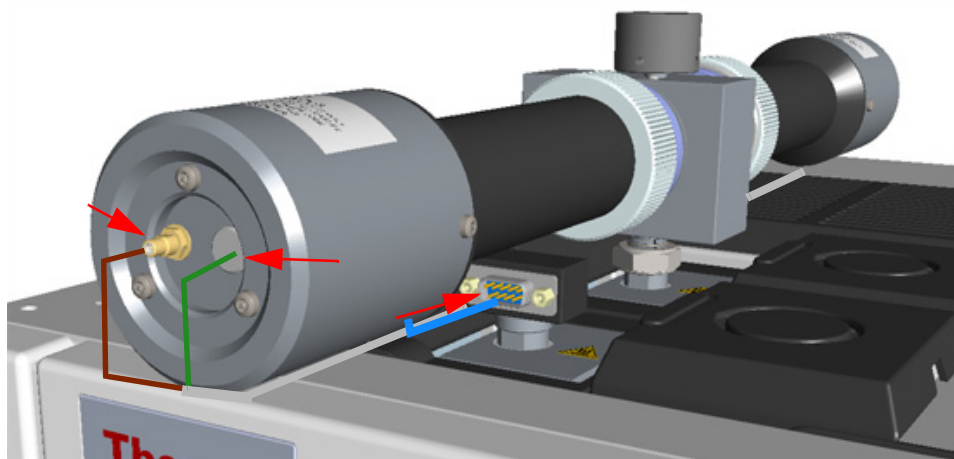

Excitation Voltage Cable   ■   Signal Cable   ■   Ignition Heating Cable   ■

#### Restart the GC and set the working parameters

16. If other external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.



17. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.

**Set the working parameters**

18. Set the parameters for the primary FPD detector (back).
  - a. Set the flow of the gases, the temperatures, the ignition threshold, and the PMT voltage for the primary photomultiplier tube according to the working analytical conditions.
19. Set the parameters for the secondary FPD detector (front).
  - a. Set the flow of the gases and the temperatures to **Off**.
  - b. Set the PMT voltage for the secondary photomultiplier tube.
20. Turn on the Flame On parameter of the primary FPD detector.



## Adding Modules

This chapter describes how to install any added injector, detector, external, module that is available for the TRACE 1300/TRACE 1310. See the *TRACE 1300 and TRACE 1310 Spare Parts Guide* for information about ordering the equipment in this chapter.


### Contents

- [Adding a SSL, SSLBKF, PTV, or PTVBKF Injector Module](#)
- [Adding a GSV Injector Module](#)
- [Adding a FID, TCD/TCD In-Series, ECD, or FPD Detector Module](#)
- [Adding a NPD Detector Module](#)
- [Adding an Aux Temperature/Cryo Module](#)
- [Adding a Helium Saver Injector Module](#)
- [Adding a PDD Module](#)
- [Adding a Generic Detector Interface](#)
- [Adding an Analog Output Interface](#)

# Adding a SSL, SSLBKF, PTV, or PTVBKF Injector Module

This section provides instructions for adding a front/back SSL, SSLBKF, PTV, or PTVBKF injector module. According to the configuration of your TRACE 1300/TRACE 1310, a dummy module is present into the free site where the injector module is not installed. See [Figure 438](#).

**Figure 438.** Add a New Front/Back Injector Module

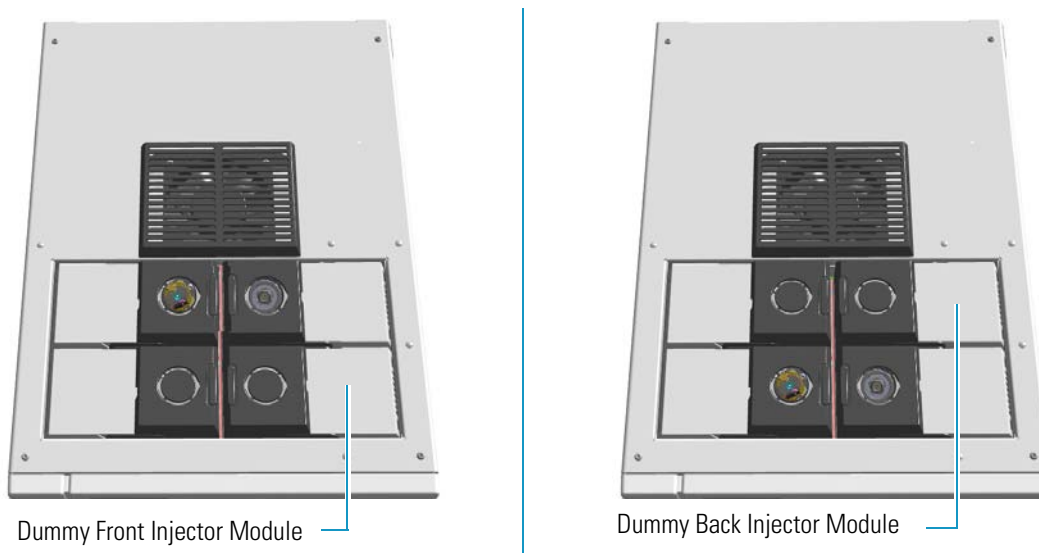

**ATTENTION** Where a dummy module is installed, the gas connection is blocked by a plug.

### ❖ To add an SSL, SSLBKF, PTV, or PTVBKF injector module


1. Put the GC in standby condition.
2. Cool the oven, injector, and detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.

6. Put the autosampler away if present.
7. Remove the dummy module from the position where the injector module will be installed.
  - a. Open the module flap cover.
  - b. Using a T20 Torxhead screwdriver, unscrew the three captive fixing screws.
  - c. Keeping the dummy module flap cover open, lift up the module from its seat into the injector housing. Place the dummy module on a clean surface.
  - d. Remove the gas block plug from the gas connections by unscrewing its fixing screw using a T20 Torxhead screwdriver. See [Figure 439](#).

**Figure 439.** Injector Gas Block Plug

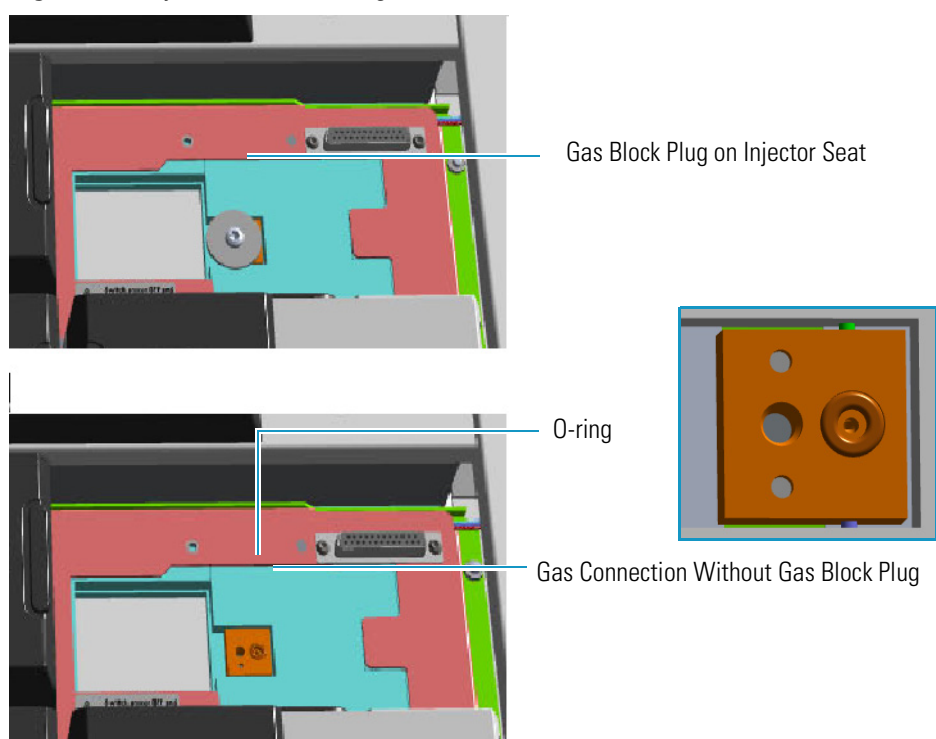

**WARNING** Make sure the o-ring is placed into its seat on the gas connection. See [Figure 439](#). **Do not install the module if the o-ring is missing.**

8. Open the front door of the GC.
9. Plug the injector module into the main frame.
  - a. Open the module flap cover.
  - b. Keeping the module flap cover open, place the module in its seat. Be sure to insert the 25-pin male connector, on the bottom of the module, into the 25-pin female connector on the injector seat of the injector housing.

## 9 Adding Modules

Adding a SSL, SSLBKF, PTV, or PTVBKF Injector Module

- c. Use a T20 Torxhead screwdriver to tighten the three captive fixing screws without overtightening.



**ATTENTION** To maintain the correct alignment the screws must be tightened in turn. Tighten each screw only a small amount before moving to the next screw. Repeat until all are secure.

- d. Close the module flap cover.
  - e. Keep the plug connected to the bottom.
10. Open the gas supplies.
  11. Check the gas supply for leaks.
    - a. Use a handheld electronic leak detector (Thermo Scientific GLD Pro leak detector or equivalent) to check each fitting for leaks.
    - b. If you detect a leak, tighten the connection and retest it.
    - c. Repeat this process until all connections are leak free.
  12. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
  13. Power on the GC.
    - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
    - b. Flip up the power switch (breaker) to the position I.
  14. Pressurize the module with the carrier gas.
  15. Check the module gas connections for leaks.
  16. Remove the plug from the bottom.
  17. Connect the analytical column end to the injector and verify the connection point.
    - a. Position the column so that the end of the column extends the proper distance above the end of the ferrule as reported in [Table 16](#).

**Table 16.** Column Insertion Depth For SSL, SSLBKF, PTV, and PTVBKF Injectors

| SSL              | SSLBKF           | PTV   | PTVBKF |
|------------------|------------------|---|--------|
| 5 mm (splitless) | 5 mm (splitless) | 30 mm   | 30 mm  |
| 10 mm (split)    | 10 mm (split)    | As far as possible into the bottom when the PTV is used as an On-Column injector. |        |


18. Close the front door of the GC.

19. If present, move the autosampler towards the module to restore the original alignment.

## Adding a GSV Injector Module

This section provides instructions for adding a front/back Gas Sampling Valve injector module. According to the configuration of your TRACE 1300/TRACE 1310, a dummy module is present into the free site where the injector module is not installed. See [Figure 440](#).

**Figure 440.** Add a New Front/Back GSV Injector Module

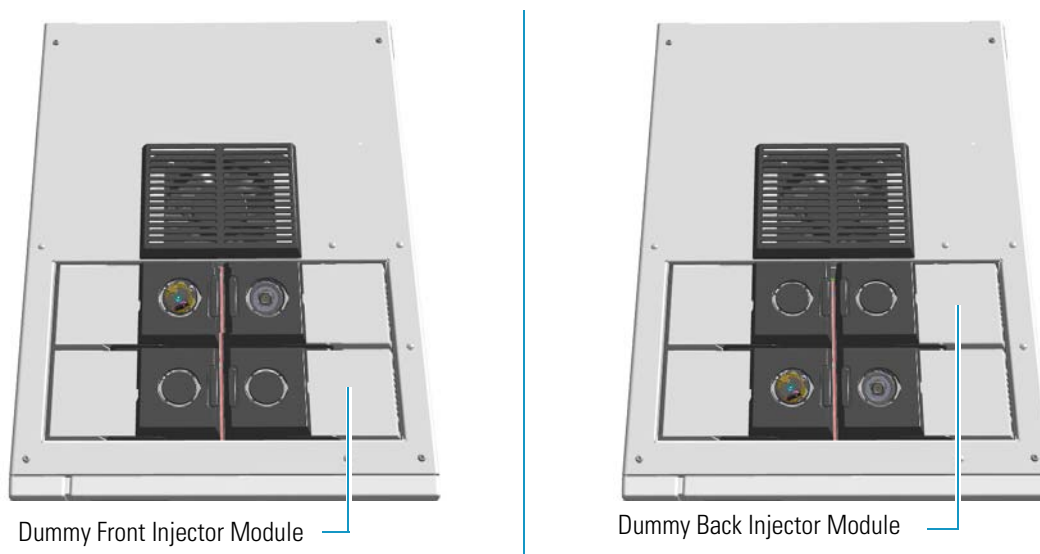

**ATTENTION** Where a dummy module is installed, the gas connection is blocked by a plug.

### ❖ To connect the gas sampling valve module to the TRACE 1300/1310 GC


1. Put the GC in standby condition.
2. Cool the oven, injector, and detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.

- b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
6. Put the autosampler away if present.
7. Remove the dummy module from the position where the injector module will be installed.
  - a. Open the module flap cover.
  - b. Using a T20 Torxhead screwdriver, unscrew the three captive fixing screws.
  - c. Keeping the dummy module flap cover open, lift up the module from its seat into the injector housing. Place the dummy module on a clean surface.
  - d. Remove the gas block plug from the gas connections by unscrewing its fixing screw using a T20 Torxhead screwdriver. See [Figure 441](#).

**Figure 441.** Injector Gas Block Plug

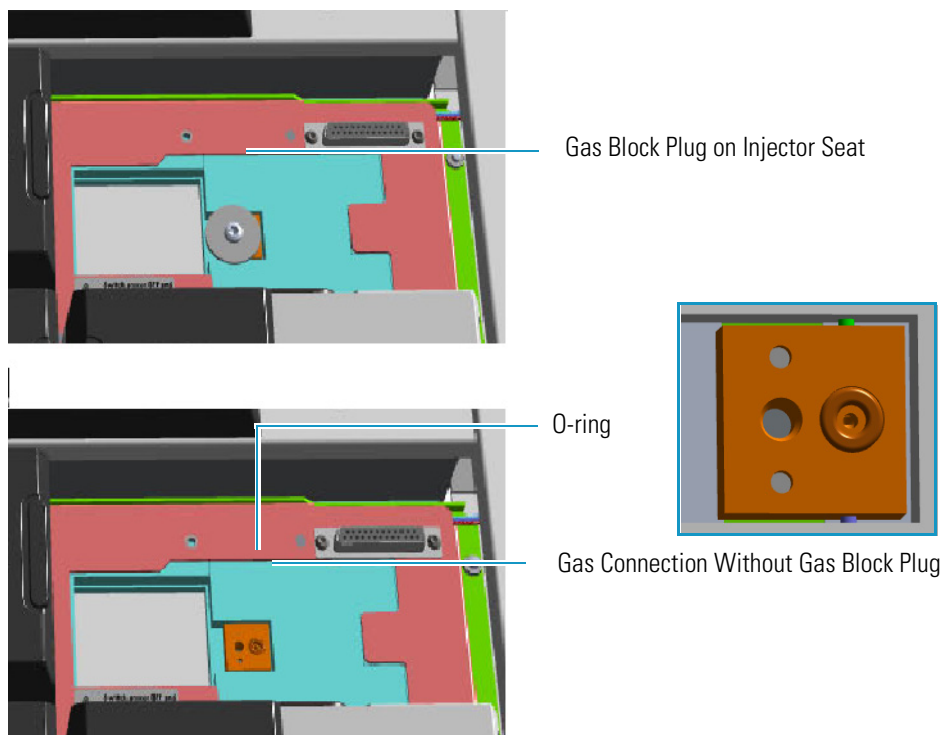

**WARNING** Make sure the o-ring is placed into its seat on the gas connection. See [Figure 439](#). **Do not install the module if the o-ring is missing.**

8. Open the front door of the GC.
9. Plug the Gas Sampling Valve module into the main frame.



- a. Open the module flap cover.
- b. Keeping the module flap cover open, place the module in its seat. Be sure to insert the 25-pin male connector, on the bottom of the module, into the 25-pin female connector on the detector seat of the detector housing.
- c. Use a T20 screwdriver to tighten the three captive fixing screws without overtightening.


**Note** Tighten the center screw first and then secure the side screws.

- d. Close the module flap cover.
10. Connect the analytical column end insert it as far as goes and withdrawn about 2 -3 mm. See the section “[Installing the Column the First Time](#)” on [page 57](#).

**Note** If you need the **backflush**, connect the **backflush system** into the GC oven following the instructions provided in the section **Connecting the GSV Backflush System** in Chapter 8 of the *TRACE 1300/TRACE 1310 Hardware Manual*.


11. Close the front door of the GC.
12. Remove the protective caps from the Sample In and Sample Out fittings. See [Figure 442](#).

**Figure 442.** Protective Caps

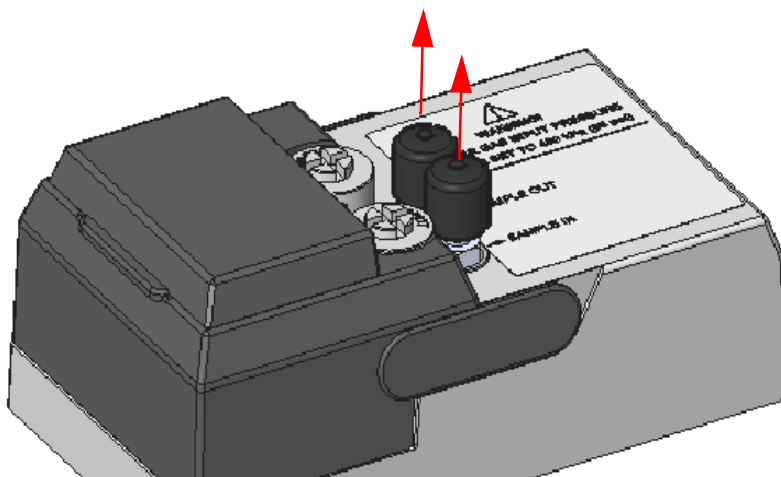


13. Connect the Sample In and Sample Out lines.
  - a. By using the proper 1/8-in. tubing, nut and ferrule, connect the inlet sample line to the **Sample IN** port on the GSV module. See [Figure 443](#).

**Figure 443.** Sample In Line Connection

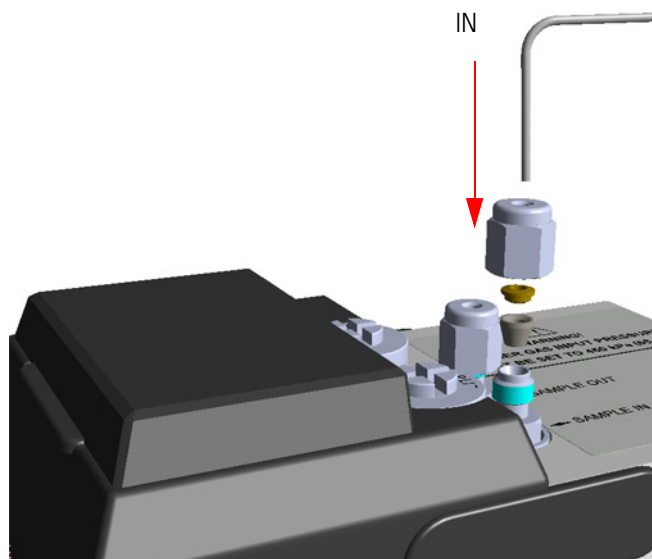

- b. By using the proper 1/8-in. tubing, nut and ferrule, connect the **Sample OUT** port on the GSV module to an exhaust system. See [Figure 444](#).


**Figure 444.** Sample Out Line Connection

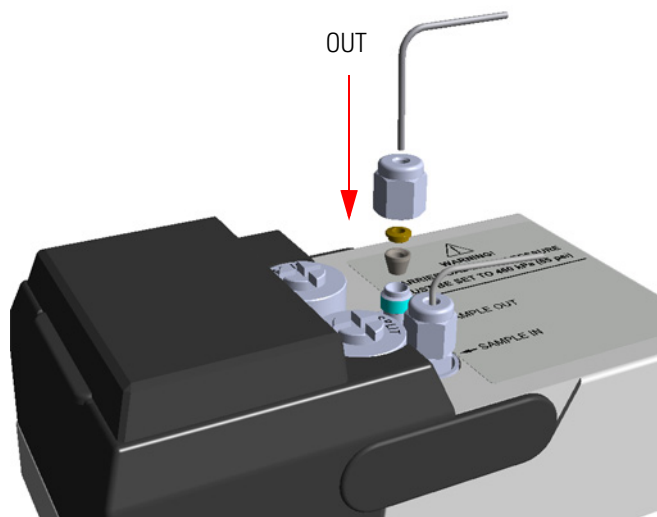

- 14. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
- 15. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
- 16. Open the module flap cover.

17. Remove the relief pins from the valve body.

**IMPORTANT** The valve is delivered with relief pins. These pins **MUST BE** removed from the valve body before working with the valve.

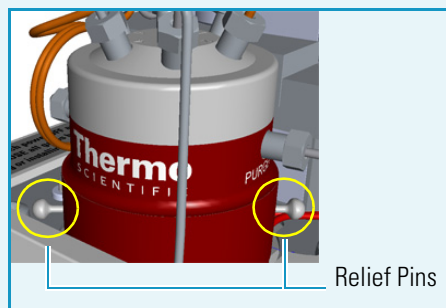

**WARNING** DO NOT SET THE CARRIER GAS INPUT PRESSURE HIGHER THAN 450 kPa (65 psig; 4.5 bar), OTHERWISE THE VALVE MAY BE DAMAGED.


**WARNING** DO NOT USE HYDROGEN AS CARRIER GAS. The module is not compatible with Hydrogen carrier gas. The same gas used as carrier gas is used to purge and to actuate the valve and must be an inert gas.

- a. Open the gas supplies.
- b. Adjust the gas pressure to 450 kPa (65 psig; 4.5 bar), actuate the valve and remove the pins.

You may need to increase the pressure a little to ease hand removal of the pins. Remember, make sure to properly readjust the actuating operating pressure after removing the pins, if you had to change it. Keep these pins in a safe place.

18. For critical applications, where air diffusion must be the lowest possible, we suggest to replace the purge vent tubing with the longer one provided. See [Figure 445](#).

**Figure 445.** Purge Vent Tubing

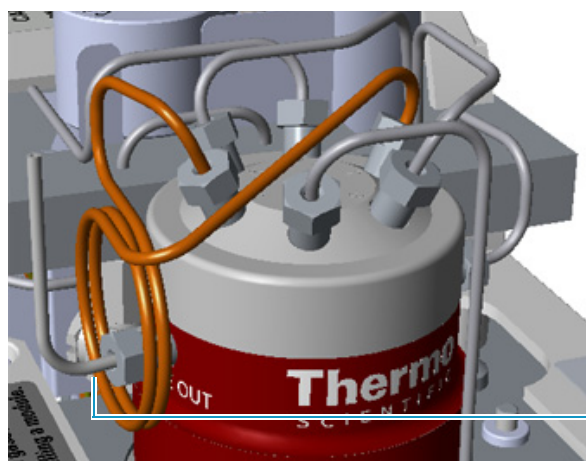

Purge Vent Tubing

## 9 Adding Modules


### Adding a GSV Injector Module

- a. Remove the purge valve tubing from the **Purge Out** port on the valve body.
  - b. Connect to the 24-in. long tubing from the Purge port on the valve body by using the proper nut and ferrule.
19. Close the module flap cover.
20. If present, move the autosampler towards the to restore the original alignment.

## Adding a FID, TCD/TCD In-Series, ECD, or FPD Detector Module

This section provides instructions for adding a front/back FID, TCD/TCD In-Series, ECD, or FPD detector module. According to the configuration of your TRACE 1300/TRACE 1310, a dummy module is present into the free site where the detector module is not installed. See [Figure 446](#).

**Figure 446.** Add a Front/Back Detector Module

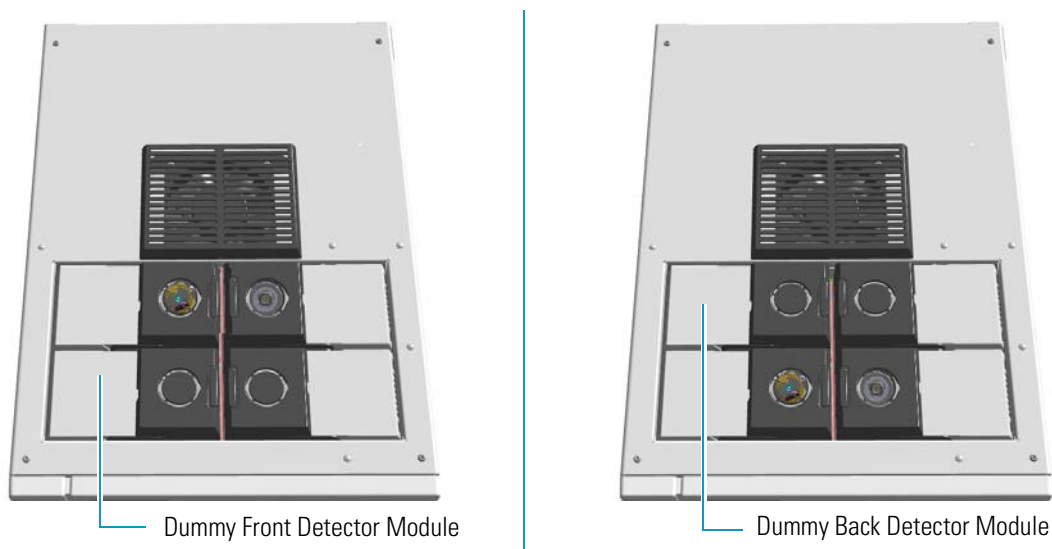

**ATTENTION** Where a dummy module is installed, the gas connection is blocked by a plug.

### ❖ To add a new FID, TCD/TCD In-Series, ECD, or FPD detector module

1. Put the GC in standby condition.
2. Cool the oven, injector, and detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.

## 9 Adding Modules

Adding a FID, TCD/TCD In-Series, ECD, or FPD Detector Module


6. Remove the dummy module from the position where the detector module will be installed.



**ATTENTION** The FPD detector must be installed in the **back** site of the detector housing. For details about the installation of the FPD detector module, see the section “[Installing the FPD Detector](#)” on [page 176](#).

- a. Open the module flap cover.
- b. Using a T20 Torxhead screwdriver, unscrew and remove the two captive fixing screws.
- c. Keeping the dummy module flap cover open, lift up the module from its seat in the injector/detector housing. Place the dummy module on a clean surface.
- d. Remove the gas block plug from the gas connections by unscrewing its fixing screw using a T20 Torxhead screwdriver. See [Figure 447](#).

**Figure 447.** Detector Gas Block Plug

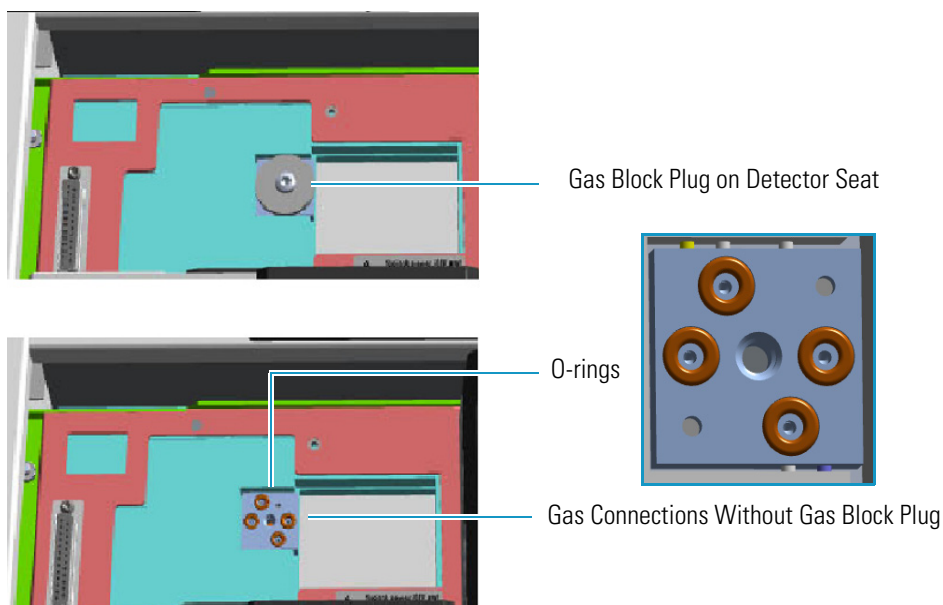

**WARNING** Make sure all the four o-rings are placed into their seats on the gas connection. See [Figure 447](#).

**Do not install the module if the o-rings are missing.**

7. Open the front door of the GC.
8. Plug the detector module into the main frame.
  - a. Open the module flap cover.

- b. Keeping the module flap cover open, place the module in its seat. Be sure to insert the 25-pin male connector, on the bottom of the module, into the 25-pin female connector on the detector seat of the detector housing.
- c. Use a T20 Torxhead screwdriver to tighten the three captive fixing screws without overtightening.



**ATTENTION** To maintain the correct alignment the screws must be tightened in turn. Tighten each screw only a small amount before moving to the next screw. Repeat until all are secure.

- d. Close the module flap cover.
9. Remove the plug and connect the analytical column end to the detector.
    - a. Position the column so that the end of the column extends the proper distance above the end of the ferrule as reported in [Table 17](#).

**Table 17.** Column Insertion Depth For FID, NPD, TCD, ECD, and FPD Detectors

| FID  | NPD | TCD | ECD   | FPD    |
|--|-----|-----|-------|--------|
| Insert the column as far as goes and withdrawn about 2 -3 mm |     |     | 23 mm | 125 mm |

**Note** If you are installing an **In-Series Connection TCD** module (TCD module modified for the connection in-series with a second detector or with a methanizer), you must connect the vent line metal tubing, that protrudes into the GC oven, to a FIDsecond detector proceeding as follows:

- Being careful not to create too narrow angles, bend the vent line metal tubing until its end reaches the bottom of the second detector.
- Making sure that the end of the vent line metal tube is perfectly straight, insert it like as an analytical column through the proper detector retaining nut and ferrule.

**ATTENTION** If you are connecting the In-Series Connection TCD module to a **FID**, **NPD**, or **FPD**, position the vent line metal tubing so that the end of the tubing extends the proper distance above the end of the ferrule as reported in [Table 17](#).

If you are connecting the In-Series Connection TCD module to an **ECD**, the insertion depth of vent line metal tubing must be **21 mm**.

10. Open the gas supplies.
11. Check the gas supply for leaks.
  - a. Use a handheld electronic leak detector (Thermo Scientific GLD Pro leak detector or equivalent) to check each fitting for leaks.
  - b. If you detect a leak, tighten the connection and retest it.

## 9 Adding Modules

Adding a FID, TCD/TCD In-Series, ECD, or FPD Detector Module

- c. Repeat this process until all connections are leak free.
12. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
13. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
14. Check the module gas connections for leaks.
15. Close the front door of the GC.




## Adding a NPD Detector Module

This section provides instructions for adding a front/back NPD detector module. According to the configuration of your TRACE 1300/TRACE 1310, a dummy module is present into the free site where the detector module is not installed. See [Figure 448](#).



**ATTENTION** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.

**Figure 448.** Add a Front/Back NPD Module

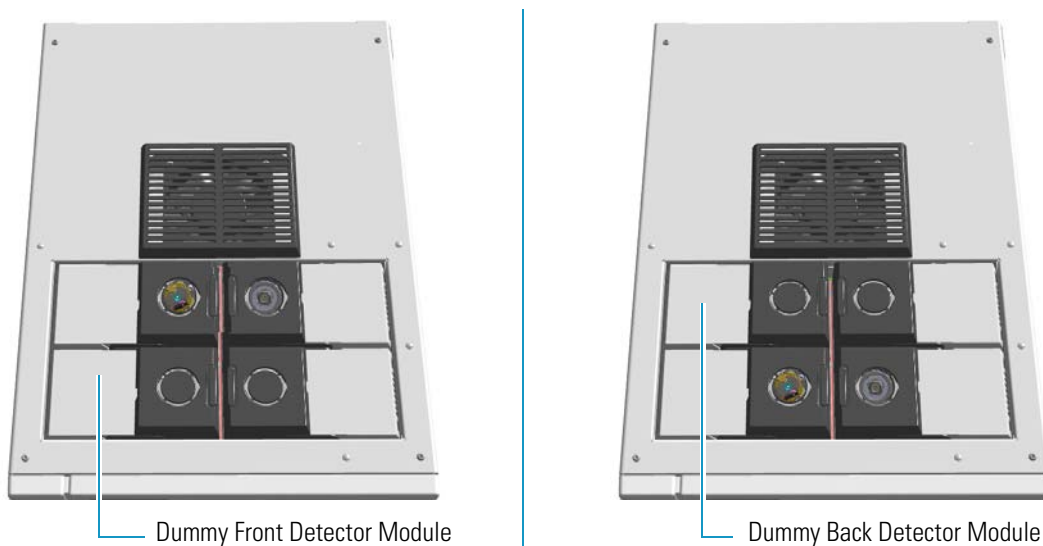

**ATTENTION** Where a dummy module is installed, the gas connection is blocked by a plug.


The NPD detector module addition also requires the installation of the NPD Thermionic Source Power Module. See [Figure 449](#).

**Figure 449.** NPD Thermionic Source Power Module



The module includes the following connections. See [Figure 450](#).

**Figure 450.** NPD Thermionic Source Power Module Connections

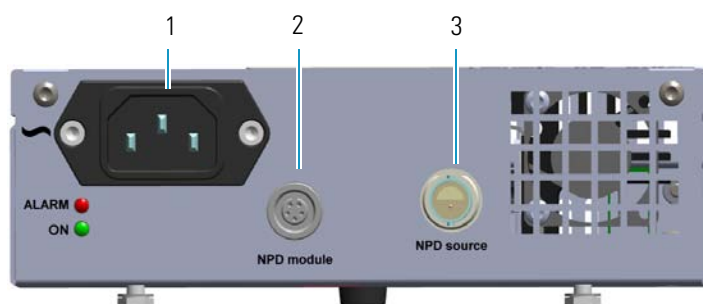

1. AC input connector.  
The module has a power rating of 120/230 Vac; 50/60 Hz; 65VA.
2. A 5-pin connector marked **NPD Module (2)** for the connection of the signal cable coming from the NPD detector module.
3. A 2-pin connector marked **NPD Source (3)** for the connection of the thermionic source assembly cable coming from the NPD detector module.

One protecting fuse is present inside the module. See [Table 18](#).

**Table 18.** Aux Temperature/Cryo Module Protecting Fuses

| Fuse | Type                  | Protections   |
|------|-----------------------|---|
| F1   | F4A 250V; (5 x 20 mm) | Circuit for the generation of the current for the thermionic source |

❖ **To add a new NPD detector module**


1. Put the GC in standby condition.
2. Cool the oven, injector, and detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
6. Remove the dummy module from the position where the detector module must be installed.
  - a. Open the module flap cover.

- b. Using a T20 Torxhead screwdriver, unscrew and remove the two captive fixing screws.
- c. Keeping the dummy module flap cover open, lift up the module from its seat in the injector/detector housing. Place the dummy module on a clean surface.
- d. Remove the gas block plug from the gas connections by unscrewing its fixing screw using a T20 Torxhead screwdriver. See [Figure 451](#).

**Figure 451.** Detector Gas Block Plug

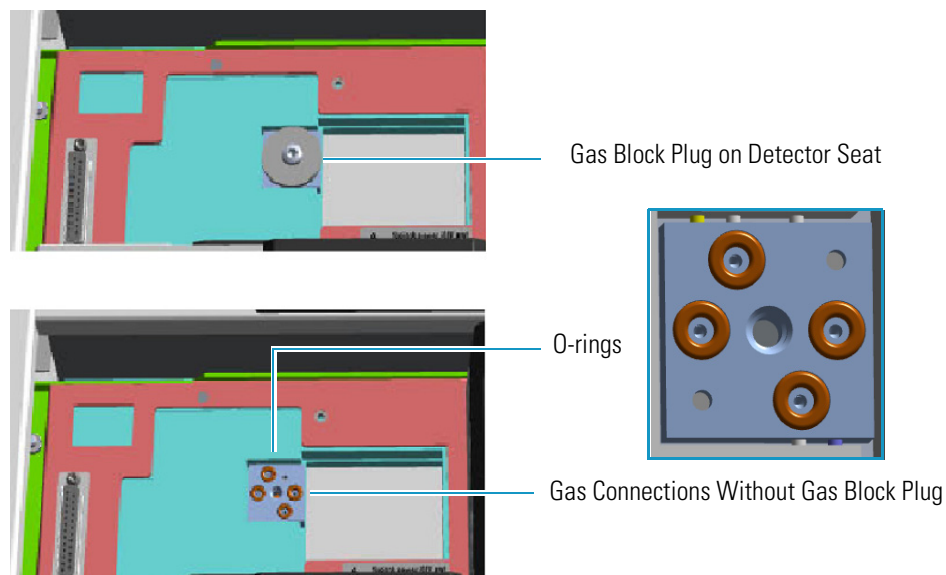

**WARNING** Make sure all the four o-rings are placed into their seats on the gas connection. See [Figure 451](#).

**Do not install the module if the o-rings are missing.**

7. Open the front door of the GC.
8. Plug the detector module into the main frame.
  - a. Open the module flap cover.
  - b. Keeping the module flap cover open, place the detector module in its seat. Be sure to insert the 25-pin male connector on the bottom of the module into the 25-pin female connector on the detector seat of the detector housing.
  - c. Use a T20 Torxhead screwdriver to tighten the three captive fixing screws without overtightening.




**ATTENTION** To maintain the correct alignment the screws must be tightened in turn. Tighten each screw only a small amount before moving to the next screw. Repeat until all are secure.

- d. Close the module flap cover.

9. Install the NPD Thermionic Source Power Module.


The module should be installed into an external module housing provided on the back of the GC. See [Figure 452](#).

**Figure 452.** External Modules Housing

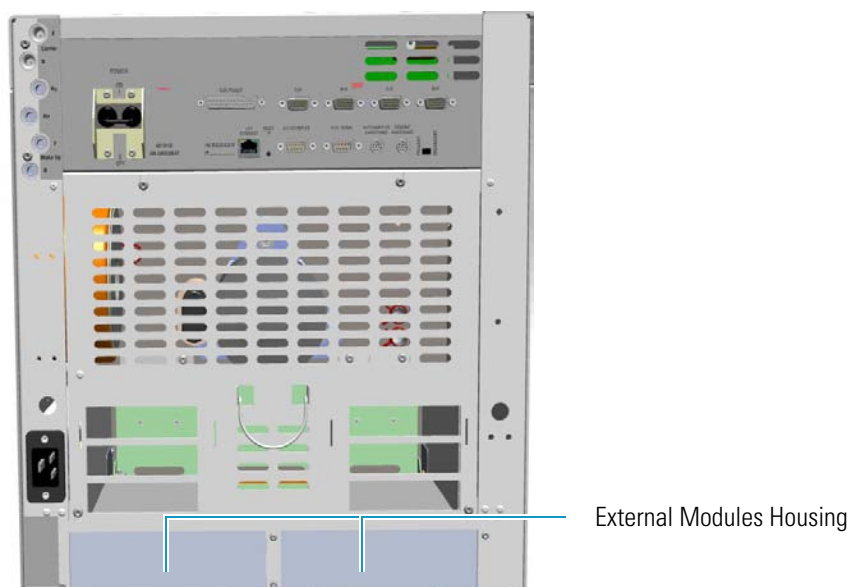

10. Remove the cover of the external modules housing where installing the module. See [Figure 453](#).

**Figure 453.** Housing Cover Removal

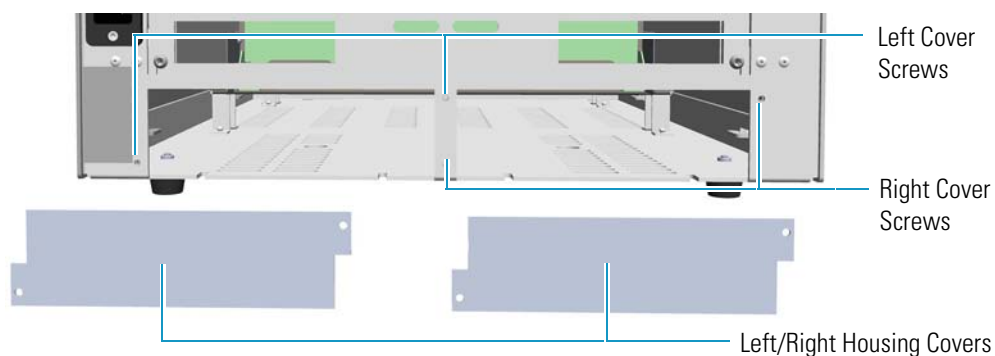


- a. Using a T20 Torxhead screwdriver, unscrew and remove the left or right housing cover screws.
  - b. Remove the covers from the housing.
11. Install the module into the housing.
    - a. Loosen the two hexagonal screws under the module. See [Figure 454](#).

**Figure 454.** Module Installation (1)



- b. Place the module into the left or right housing until the hexagonal screws couple the slots on the floor of the GC. See [Figure 455](#).

**Figure 455.** Module Installation (2)

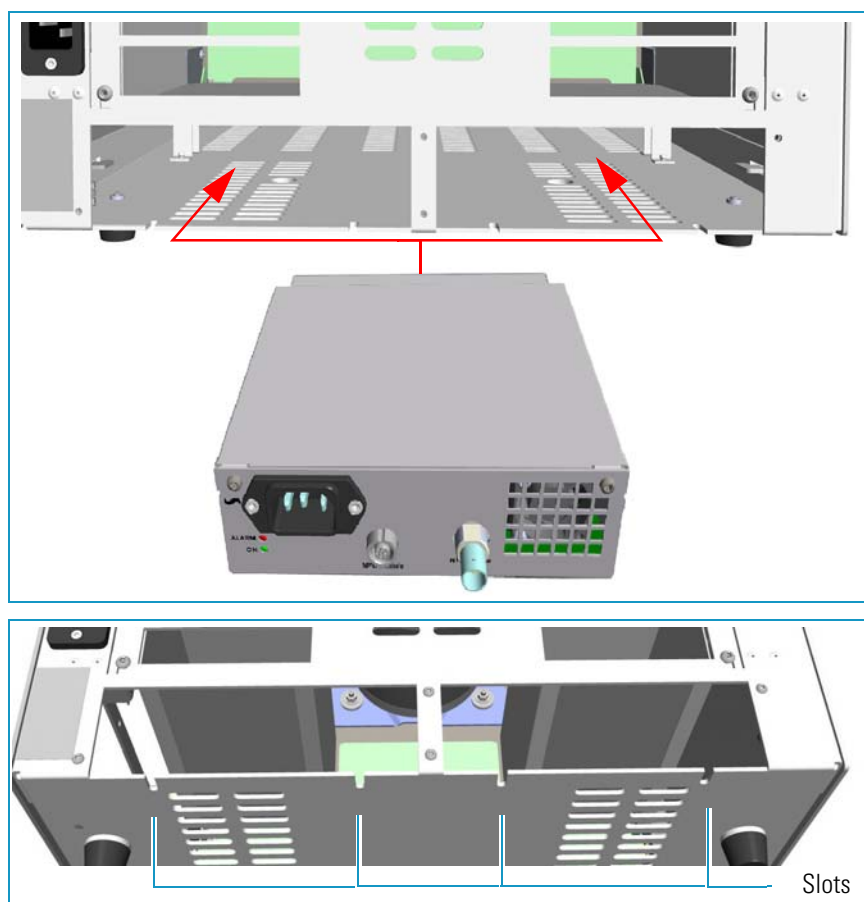


- c. Fingertighten the hexagonal screws slightly, or using a 10-mm wrench.

**Note** Always keep the hexagonal screws in their place. This allows you an easier removal of the auxiliary module when necessary.

The result of the installation is shown in [Figure 456](#).

**Figure 456.** NPD Thermionic Source Power Module Installed into the GC

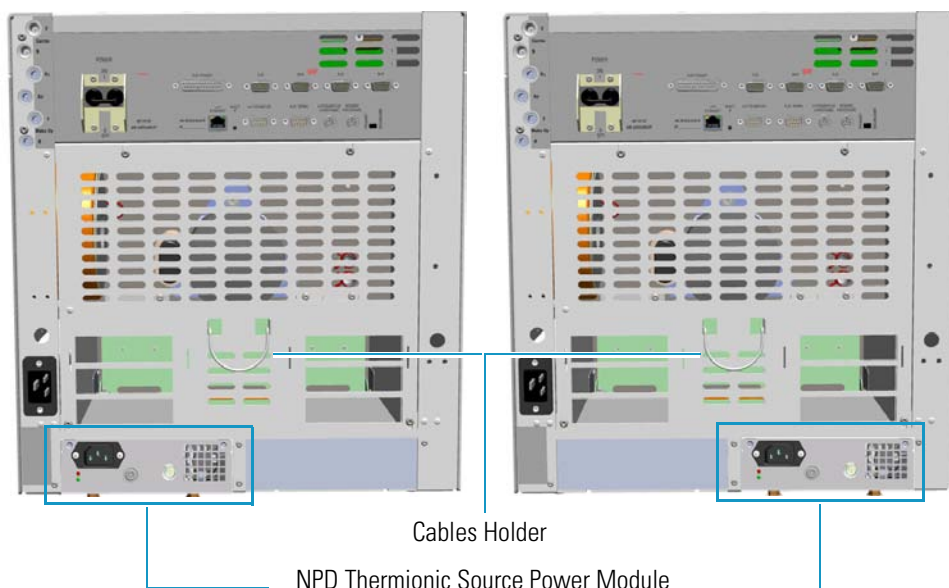


12. Connect the NPD Thermionic Source Power Module.

**Tip** To avoid the contact with the hot air from the vents, it is suggested to gather the electrical cables into cables holder.


- Open the module flap cover.
- Using the cable provided, connect the source assembly connector on the NPD Thermionic Source Power Module to the thermionic source and twist the ring to lock the connection. See [Figure 457](#).

**Figure 457.** Thermionic Source Assembly Cable Connection

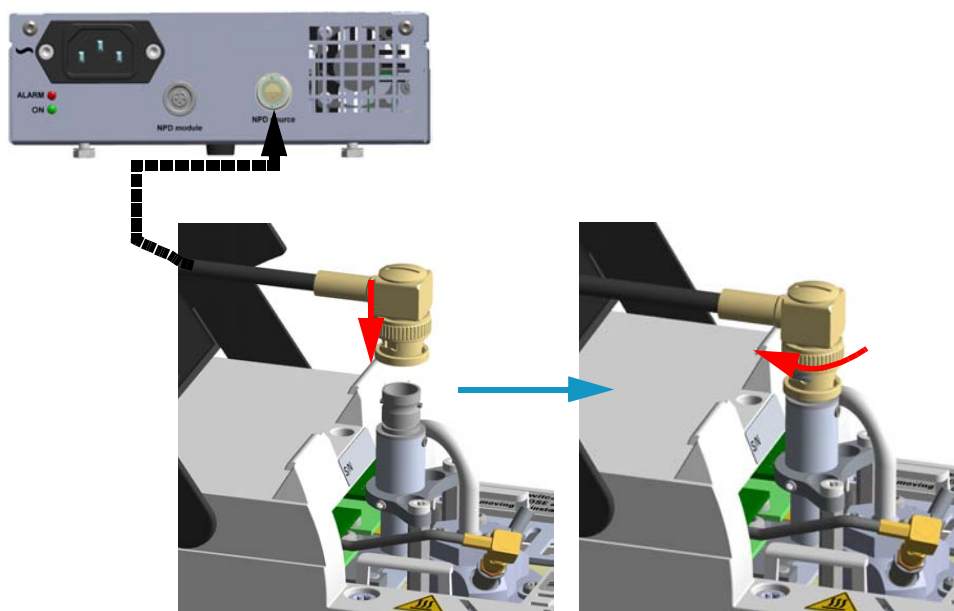

- c. Using the cable provided, connect the signal connector on the NPD Thermionic Source Power Module to the connector on the detector module. See [Figure 458](#).

**Figure 458.** Signal Cable Connection

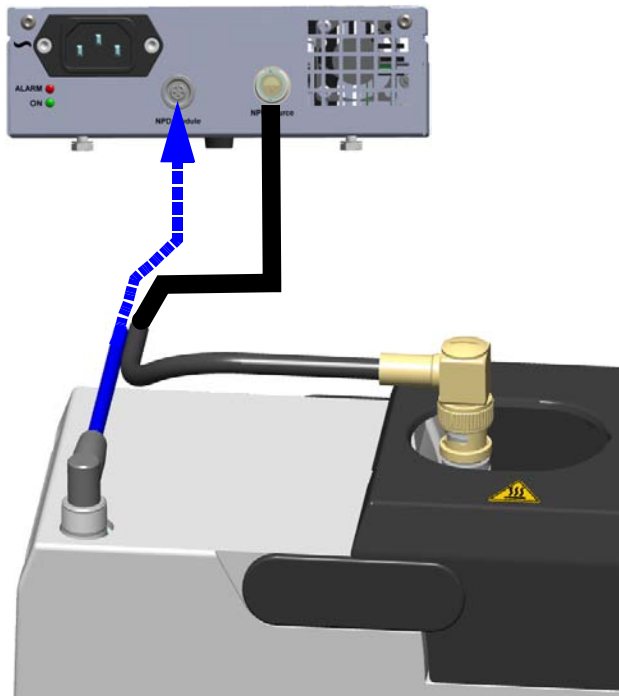

13. Remove the plug and connect the column end to the detector.
14. Open the gas supplies.
15. Check the gas supply for leaks.
  - a. Use a handheld electronic leak detector (Thermo Scientific GLD Pro leak detector or equivalent) to check each fitting for leaks.
  - b. If you detect a leak, tighten the connection and retest it.
  - c. Repeat this process until all connections are leak free.
16. Plug the power cable to the AC input connector of the NPD Thermionic Source Power Module, and to the wall outlet. The LED marked **On** blinks green.



**IMPORTANT** The NPD Thermionic Source Power Module must be powered at the same line voltage of the main GC system.

If other external modules are present, plug the power cable to the AC Input connector of each module, and to the wall outlet.

17. Power on the GC.

## 9 Adding Modules

### Adding a NPD Detector Module

- a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
18. Check the module gas connections for leaks.
19. Close the front door of the GC.
20. Set the make-up gas on.
21. Switch on the thermionic source. Hydrogen and air are automatically opened.  
The LED marked **On** of the NPD Thermionic Source Power Module becomes solid green.
22. With all gases on, heat the detector to 150 °C and hold for about 15 minutes, then increase the temperature to 250 °C and hold for 15 minutes.
23. Increase the temperature to operating value (300 °C to 320 °C recommended). Allow 15 minutes for equilibration.
24. Check the NPD leakage current. If > 2.0 pA, verify the installation.



## Adding an Aux Temperature/Cryo Module

This section provides instructions for updating your TRACE 1300/TRACE 1310 with the Aux Temperature/Cryo module. See [Figure 459](#).




**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.

**Figure 459.** Aux Temperature/Cryo Module



The module includes the following connections. See [Figure 460](#).

**Figure 460.** Aux Temperature/Cryo Module Connections

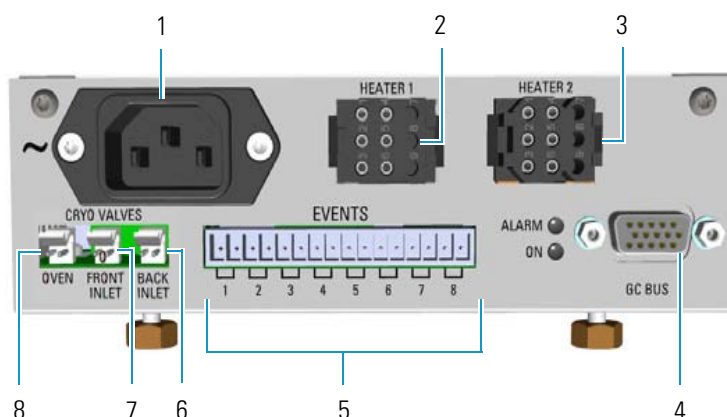

1. AC input connector.  
The module has a power rating of 120/230 Vac; 50/60 Hz; 355 VA.
2. A 6-pin connector marked **Heater 1** for the connection of the transfer line of a mass spectrometer.
3. A 6-pin connector marked **Heater 2** for the connection of the transfer line of a second mass spectrometer.
4. 15-pin female connectors marked **Bus** for the communication with the GC.
5. A 16-pin connector marked **Events** for the connection of **eight** external solenoid valves (24 V, 2 W max). The valves can be programmed individually during the running of the instrument.

## 9 Adding Modules

### Adding an Aux Temperature/Cryo Module

6. A 2-pin connector marked **Back Inlet** for the connection of the solenoid valves for the back PTV/PTVBKF cryogenic system.
7. A 2-pin connector marked **Front Inlet** for the connection of the solenoid valves for the front PTV/PTVBKF cryogenic system.
8. A 2-pin connector marked **Oven** for the connection of the solenoid valves for the Oven cryogenic system.

Four protecting fuses are present inside the module. See [Table 19](#).

**Table 19.** Aux Temperature/Cryo Module Protecting Fuses


| Fuse | Type                    | Protections          |
|------|-------------------------|----------------------|
| F1   | F2A 250V; (5 x 20 mm)   | Heater 1             |
| F2   | F2A 250V; (5 x 20 mm)   | Heater 2             |
| F3   | T2A 250V; (5 x 20 mm)   | 24 V                 |
| F4   | F1.6A 250V; (5 x 20 mm) | 24 V solenoid valves |

The module features the following:

- The valve control of the Oven cryo option.
- The valve control of the front/back PTV cryo option.
- Up two heater controls for the DSQ II, ITQ, or TSQ Quantum transfer line.
- Sixty-three timed events to automatically drive up to eight external on/off solenoid valves.

The module should be installed into a free external module housing provided on the back of the GC. See [Figure 461](#).

**Figure 461.** External Modules Housing

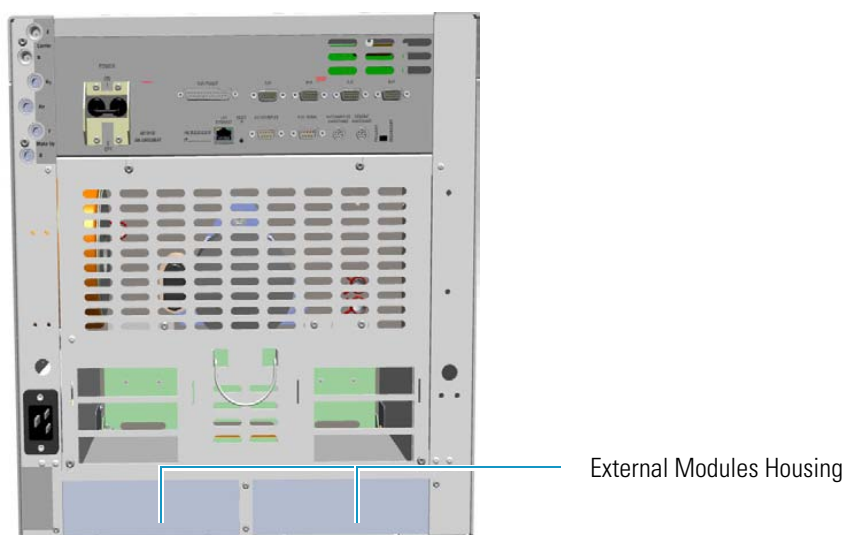


❖ **To add an Aux Temperature/Cryo Module**

1. Put the GC in standby condition.
2. Cool the oven, injectors, and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
6. Remove the cover of the external modules housing where installing the module.  
See [Figure 462](#).

**Figure 462.** Housing Cover Removal

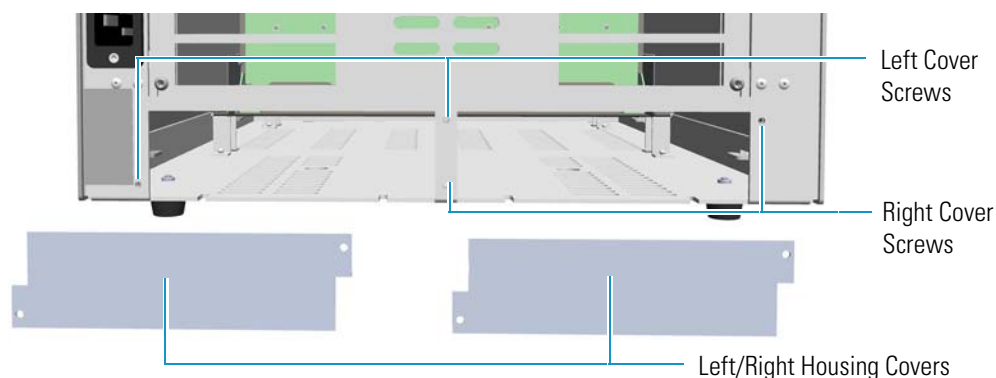


- a. Using a T20 Torxhead screwdriver, unscrew and remove the left or right housing cover screws.
  - b. Remove the covers from the housing.
7. Install the module into the housing
    - a. Loosen the two hexagonal screws under the module. See [Figure 463](#).

## 9 Adding Modules


Adding an Aux Temperature/Cryo Module

**Figure 463.** Module Installation (1)

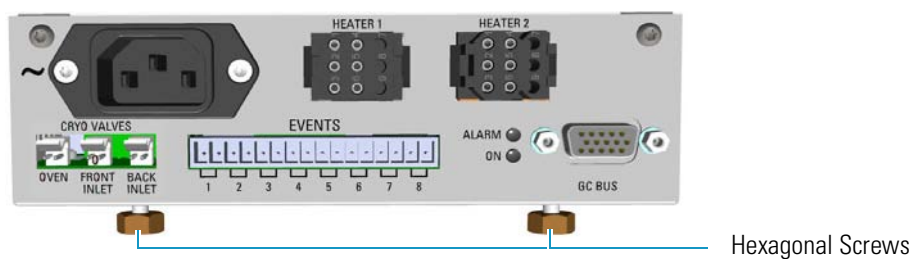

- b. Place the module into the left or right housing until the hexagonal screws couple the slots on the floor of the GC. See [Figure 464](#).

**Figure 464.** Module Installation (2)

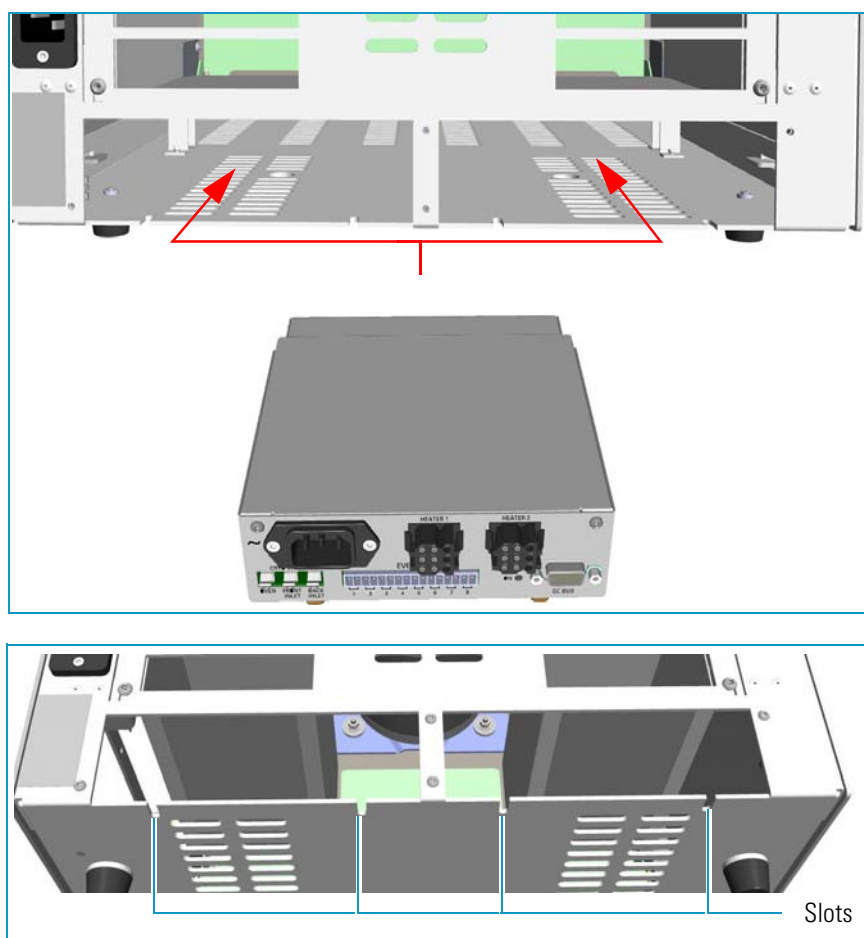


- c. Finger-tighten the hexagonal screws slightly, or using a 10-mm wrench.

**Note** Always keep the hexagonal screws in their place. This allows you an easier removal of the auxiliary module when necessary.

The result of the installation is shown in [Figure 465](#).

**Figure 465.** Aux Temperature/Cryo Module Installed into the GC

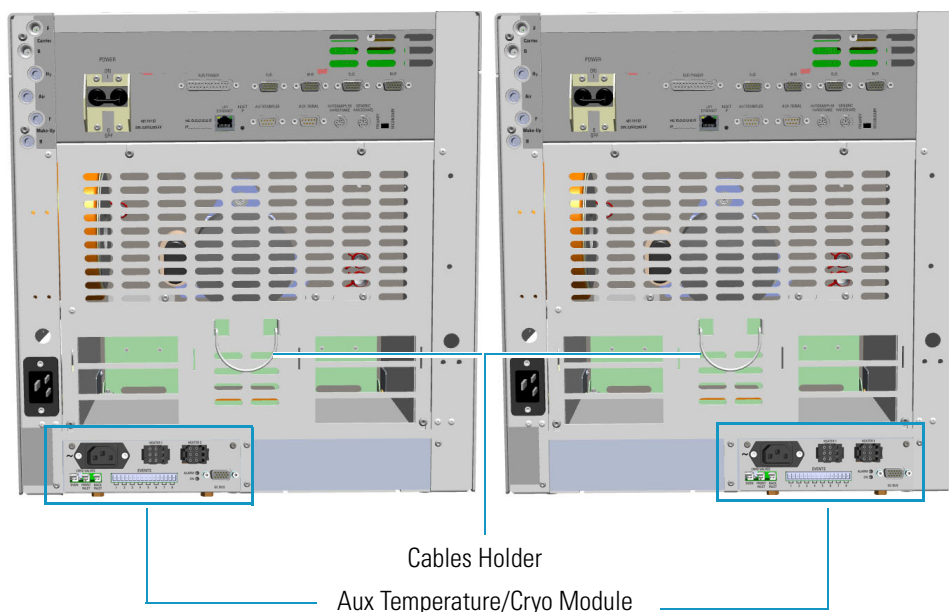

8. Connect the Aux Temperature/Cryo module.

**Tip** To avoid the contact with the hot air from the vents, it is suggested to gather the electrical cables into cables holder.

- a. Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.
  - If a DSQ II, ITQ, or TSQ Quantum mass spectrometer is coupled with the GC, connect the heater cable coming from the transfer line to the connector marked **Heater 1** or **Heater 2**.



**ATTENTION** An extension cable transfer line is needed when connecting an ITQ mass spectrometer to the Aux Temperature/Cryo module.

- If the cryo option for the Oven is present, connect the 2-pin connector marked **Cryo Valves - Oven** to the cryo solenoid valve using the cable provided.
- If the cryo option for the front/back PTV/PTVBKF injector is present, connect the 2-pin connector marked **Cryo Valves - Front Inlet** or **Cryo Valves-Back Inlet** to the cryo solenoid valve using the cable provided.
- If the on/off activation of external solenoid valves (up to eight) is required, connect a 2-pin connector marked **Event 1÷8** to each external solenoid valve.
- Plug the power cable to the AC Input connector on the front of the module and to the wall outlet. The LED marked **On** lights after the GC is powered on.



**IMPORTANT** The Auxiliary Temperature/Cryo module must be powered at the same power supply voltage required for the GC.

9. Open the gas supplies.
10. If other external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
11. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker), located at the back of the instrument, to the position I.
12. Configure and enable the system installed through the user interface of your GC, or through the CDS in use. Refer to the *TRACE 1300 and TRACE 1310 User Guide*.

## Adding a Helium Saver Injector Module

The Instant Connect Helium Saver Injector Module (HeS-S/SL) is designed as a “double wide” module that fits into the space provided by one detector module and one injector module on the upper deck of your TRACE 1300/TRACE 1310 GC.

### Installing the Module

The module may be placed in either the Front or Back position on the upper deck.

- If the module is located in the **Front** position, plumb the helium supply to the Front Make Up gas connection on the back panel of the GC.
- If the module is located in the **Back** position, plumb the helium to the Back Make Up gas connection on the back panel of the GC.

Plumb the nitrogen gas to the corresponding Front or Back Carrier gas connections.



**IMPORTANT** Although nitrogen will not be the actual carrier gas for the analytical separation, it will be the carrier gas during the injection, and it is necessary to plumb the nitrogen to the carrier gas input.

**Note** For optimal results, use a high capacity oxygen trap on the nitrogen supply, and a heated zirconium alloy gettering trap on the helium supply. These traps can be purchased together as Thermo Scientific part number 1R120577-0001.

Alternatively, you may opt to provide own helium purification. Due to the low flow rate of helium employed, conventional chemical traps (non-heated traps) may actually contaminate the gas supply. Heated zirconium-based traps specifically designed for helium are ideal. These traps can also remove nitrogen, which is difficult to eliminate from conventional traps.

❖ **To install the Instant Connect Helium Saver Injector Module**

1. Cool the GC and MS heated zones and then shut off all carrier and detector gases on the local GC user interface as well as at the source cylinder.
2. Remove the detector and injector “dummy modules” from the locations where the Helium Saver S/SL module will be installed by loosening the two captive T20 screws on each module. See [Figure 466](#).

**Figure 466.** Removing the Dummy Modules




3. Remove the blanking washers from the detector and injector pneumatic network. See [Figure 467](#).



## 9 Adding Modules

### Adding a Helium Saver Injector Module

**Figure 467.** Removing the Blanking Washers

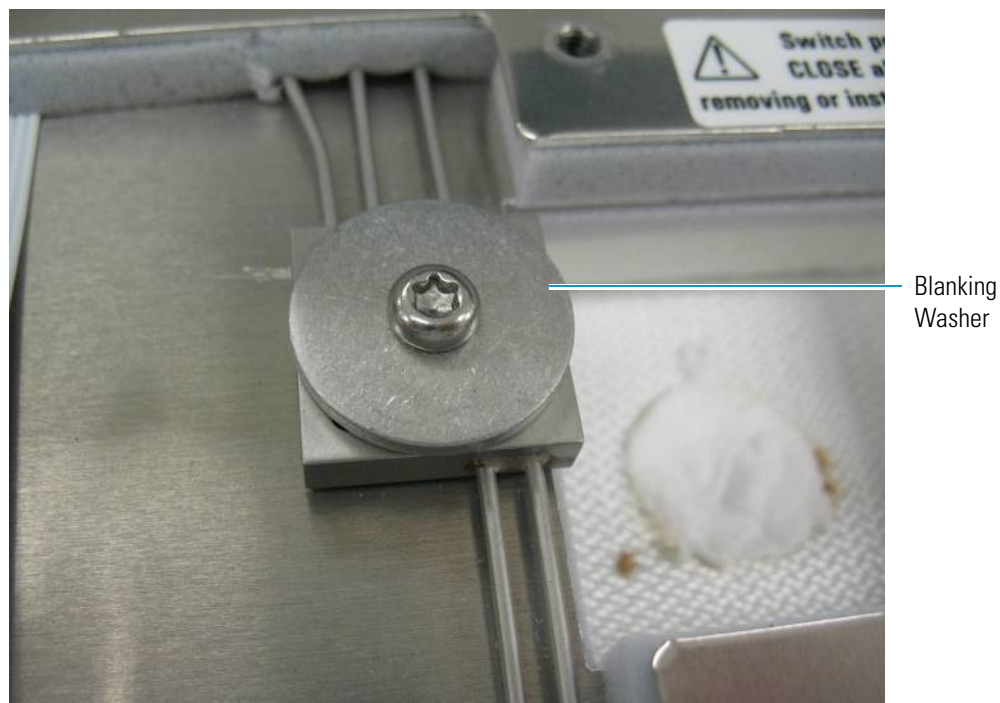

4. Ensure there is an o-ring present in the o-ring seat of each position on the pneumatic network block. See [Figure 468](#).

**Figure 468.** Confirming Presence of O-rings



5. Plumb a cylinder of high purity helium 99.999% to the appropriate make-up gas connection on the back of the GC using the adapter provided. See [Figure 469](#). The Thermo Scientific gas purification kit part number 1R120577-0001 contains a heated gettering trap which can be used directly with this fitting.



**Figure 469.** Plumbing the Helium Cylinder




Alternatively, a customer supplied gettering trap can be used along with the provided 1/16 in. stainless steel tubing. Set the regulator pressure such that helium purges the make-up line of the pneumatic network for 15 min at a flow of 20–100 mL/min.

The flow can be measured directly from the detector block at the left most port as shown in [Figure 468](#) on [page 450](#). This will purge the gas line as well as eliminate air from the gettering material of an in-line heated helium purifier. Do this before applying power to the purifier. Ensure the hose leading to the flowmeter is clean and dust free before holding it against the o-ring seat.

6. Carefully insert the Instant Connect Helium Saver Injector Module into position on the GC top deck after first inserting the ceramic insulator shown in [Figure 470](#).

**Figure 470.** Positioning the Instant Connect Helium Saver Injector Module

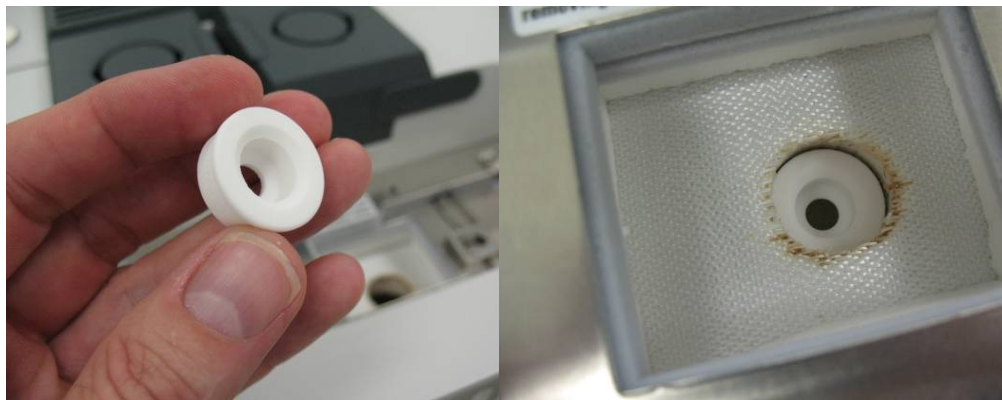

Also, be sure to remove the column nut from the injector if it has not already been removed. This prevents damage to the ceramic insulator when the module is inserted.

**Note** It may be helpful to insert the left side of the module just prior to the right side in order to clear the space properly. The module fits snugly into position but should not require undue force to clear the sheet metal opening.

## 9 Adding Modules

### Adding a Helium Saver Injector Module

- Secure the T20 screws near the center line of the module leading to the injector and detector pneumatic network prior to tightening the captive screws which hold the module. See [Figure 471](#).

**Figure 471.** Securing the Instant Connect Helium Saver Injector Module in Place



It will be necessary to push down on the right side of the module near the back edge to ensure engagement of the 25 pin D-sub electrical connector. Wiggle the module as necessary to engage the pins and then secure the screws.

- At this point, the regulator that delivers the helium can be set to the appropriate pressure: typically 110 psig (760 kPa; 7.6 bar) will suffice for 0.25 mm i.d. columns of 30 m length.

Dial the pressure up slowly while tapping the pressure gauge in order to set the appropriate pressure. [Table 20](#) serves as a guide for setting the correct helium regulator pressure.

**Table 20.** Determining Correct Helium Regulator Pressure for Column Type (Sheet 1 of 2)

| Column Length | Column i.d. | Desired Flow * | Required helium regulator pressure psig (kPa; bar) ** |
|---------------|-------------|----------------|---|
| 5 m/10 m      | 0.10 mm     | 0.4 mL/min     | 110 psig (760 kPa; 7.6 bar)                           |
|               |             | 0.5 mL/min     | 120 psig (830 kPa; 8.3 bar)                           |
|               |             | 0.6 mL/min     | 130 psig (900 kPa; 9 bar)                             |
| 10 m/20 m     | 0.18 mm     | 0.8 mL/min     | 100 psig (690 kPa; 6.9 bar)                           |
|               |             | 1.0 mL/min     | 110 psig (760 kPa; 7.6 bar)                           |
|               |             | 1.2 mL/min     | 120 psig (830 kPa; 8.3 bar)                           |
|               |             | 1.5 mL/min     | 130 psig (900 kPa; 9 bar)                             |
| 15 m          | 0.25 mm     | 2.0 mL/min     | 140 psig (970 kPa; 9.7 bar)                           |
|               |             | 1.0 mL/min     | 100 psig (690 kPa; 6.9 bar)                           |
|               |             | 1.2 mL/min     | 110 psig (760 kPa; 7.6 bar)                           |
|               |             | 1.7 mL/min     | 120 psig (830 kPa; 8.3 bar)                           |

**Table 20.** Determining Correct Helium Regulator Pressure for Column Type (Sheet 2 of 2)

| Column Length | Column i.d. | Desired Flow * | Required helium regulator pressure psig (kPa; bar) ** |
|---------------|-------------|----------------|---|
| 30 m          | 0.25 mm     | 2.2 mL/min     | 130 psig (900 kPa; 9 bar)                             |
|               |             | 2.7 mL/min     | 140 psig (970 kPa; 9.7 bar)                           |
|               |             | 1.0 mL/min     | 100 psig (690 kPa; 6.9 bar)                           |
|               |             | 1.2 mL/min     | 110 psig (760 kPa; 7.6 bar)                           |
|               |             | 1.7 mL/min     | 120 psig (830 kPa; 8.3 bar)                           |
|               |             | 2.2 mL/min     | 130 psig (900 kPa; 9 bar)                             |
| 60 m          | 0.25 mm     | 2.7 mL/min     | 140 psig (970 kPa; 9.7 bar)                           |
|               |             | 3.2 mL/min     | 150 psig (1030 kPa; 10.3 bar)                         |
|               |             | 1.0 mL/min     | 120 psig (830 kPa; 8.3 bar)                           |
|               |             | 1.2 mL/min     | 130 psig (900 kPa; 9 bar)                             |
|               |             | 1.7 mL/min     | 140 psig (970 kPa; 9.7 bar)                           |
|               |             | 2.2 mL/min     | 150 psig (1030 kPa; 10.3 bar)                         |
| 100 m         | 0.25 mm     | 1.0 mL/min     | 130 psig (900 kPa; 9 bar)                             |
|               |             | 1.2 mL/min     | 140 psig (970 kPa; 9.7 bar)                           |
|               |             | 1.7 mL/min     | 150 psig (1030 kPa; 10.3 bar)                         |
| 30 m          | 0.32 mm     | 1.5 mL/min     | 100 psig (690 kPa; 6.9 bar)                           |
|               |             | 2.0 mL/min     | 110 psig (760 kPa; 7.6 bar)                           |
|               |             | 2.5 mL/min     | 120 psig (830 kPa; 8.3 bar)                           |
|               |             | 3.0 mL/min     | 130 psig (900 kPa; 9 bar)                             |
|               |             | 3.5 mL/min     | 140 psig (970 kPa; 9.7 bar)                           |
|               |             | 1.5 mL/min     | 110 psig (760 kPa; 7.6 bar)                           |
| 60 m          | 0.32 mm     | 2.0 mL/min     | 120 psig (830 kPa; 8.3 bar)                           |
|               |             | 2.5 mL/min     | 130 psig (900 kPa; 9 bar)                             |
|               |             | 3.0 mL/min     | 140 psig (970 kPa; 9.7 bar)                           |
|               |             | 3.5 mL/min     | 150 psig (1030 kPa; 10.3 bar)                         |
|               |             | 1.5 mL/min     | 130 psig (900 kPa; 9 bar)                             |
|               |             | 2.0 mL/min     | 140 psig (970 kPa; 9.7 bar)                           |
| 100 m         | 0.32 mm     | 2.5 mL/min     | 150 psig (1030 kPa; 10.3 bar)                         |

\* For flow rates not specifically listed, round up in pressure to the next highest value. For example, if 1.5 mL/min is desired using a 0.25 mm ID column of 15 m length, use a helium regulator pressure of 110 psig (760 kPa; 7.6 bar).

\*\* Minimum pressure required to avoid nitrogen back diffusion into the helium carrier gas. Higher pressure will result in a slightly higher (than minimum) consumption of helium, but will not result in adverse analytical performance.

9. Plumb a nitrogen supply to the appropriate Carrier input at the back of the GC as shown in [Figure 472](#).

**Figure 472.** Plumbing the Nitrogen Supply to the GC



10. The Thermo Scientific gas purification kit contains a high capacity oxygen trap that should be placed in the nitrogen line. Be sure to purge the regulator and gas line before attaching the oxygen filter.
11. Keep nitrogen flowing through the trap while attaching the filter to the back of the GC in order to prevent oxygen from entering the trap.
12. Adjust the pressure regulator to supply a pressure sufficient for the analytical method in use. In general, 75 psig (517 kPa; 5.17 bar) is more than sufficient for most applications. Very small bore capillaries may require higher pressures.


## Installing the Column

The Instant Connect Helium Saver Injector Module is designed to be used with Silflow™ metallic ferrules. This allows for very accurate trimming of the column following compression of the ferrule onto the column. It is important that only 5 mm of column protrudes past the tip of the ferrule for proper operation. The inlet has been designed to work optimally with 0.25 mm i.d. columns, although larger (up to 0.32 mm i.d.) or narrower i.d.'s can be used successfully.

### ❖ To install the column

1. Insert the column through the SilFlow™ nut and ferrule as shown in [Figure 473](#).

**Figure 473.** Installing the Column Nut and Ferrule

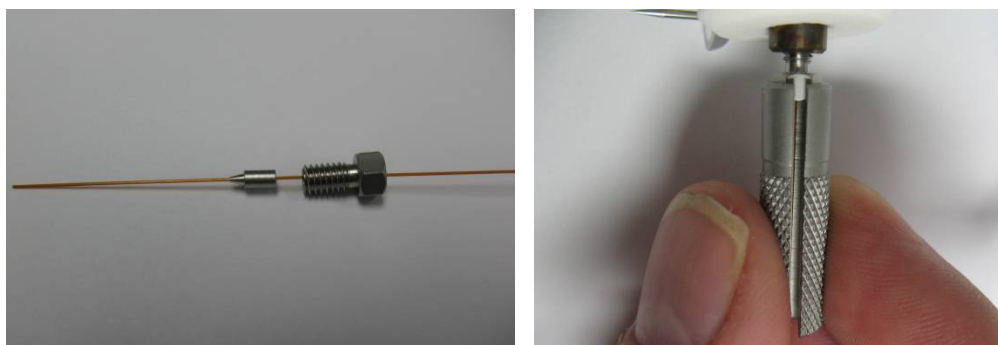

2. Allow a few centimeters of column to extend past the tip of the ferrule and insert it into the base of the inlet. It will be necessary to gently poke around in order to find the small bore cone that serves as a column guide and ferrule seat.
3. Use the knurled tool to tighten the nut by finger force only, until the ferrule grabs the column, and the column no longer slides in the bore of the ferrule.

The column connection should appear as shown in [Figure 474](#).


**Figure 474.** Installing the Column in the Injector



The column should be fully captured by the ferrule without the ability to slide it up or down.

4. At this point, remove the nut and ferrule assembly and confirm the column will not slide in the ferrule (See first **Tip** on [page 455](#)). Trim the column such that only 5 mm extends past the tip of the ferrule.
5. Carefully reinstall the column and again tighten the ferrule. The column is now installed in the injector.

**Tip** An indispensable tool for removal of the column from the inlet is a dental pick or thumbtack. When the SilFlow ferrule engages the inlet on tightening, it is slightly deformed at the tip in order for the sealing to occur. This causes the ferrule to become “stuck,” which is a normal occurrence. The ferrule can easily be removed by inserting the pointed tip of the dental pick gently along the side of the ferrule and pressing vertically to cause the ferrule to be displaced to the side. Gently pulling on the column at the same time will dislodge the ferrule.

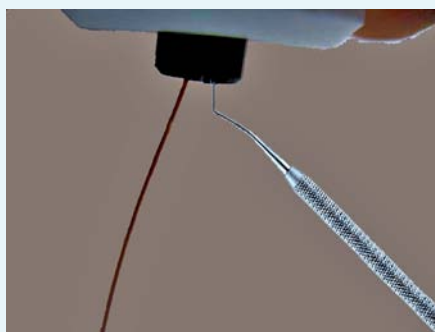

## 9 Adding Modules

### Adding a Helium Saver Injector Module

**Note** The standard outfit kit for the Instant Connect Helium Saver Injector Module is equipped with SilFlow ferrules having internal diameters of 0.50, 0.40, and 0.35 mm.

If you are using a 0.32 mm id column, use the 0.50 mm i.d. ferrules supplied.

If you are using a 0.25 mm id column or smaller, use the 0.35 mm id ferrule if the column will pass through it.


This will allow easy finger tightening of the ferrule using the hand tool without undue torque. In other words, use the smallest ferrule that will fit your columns.

The ideal situation is when the column passes through the hole with slight pressure.

The replacement Thermo part numbers are as follows:

- 0.50 mm ferrule **10ea.** 29063467 (for 0.32 columns)
- 0.40 mm ferrule **10ea.** 29063466 (for 0.25 columns and smaller on the high end of the o.d. tolerance)
- 0.35 mm ferrule **10ea.** 29063465 (for 0.25 columns and smaller on the low end of the o.d. tolerance).

**Tip** When installing the column, nut and ferrule, lift up slightly on the ceramic insulator surrounding the base of the injector insert as shown in figure.

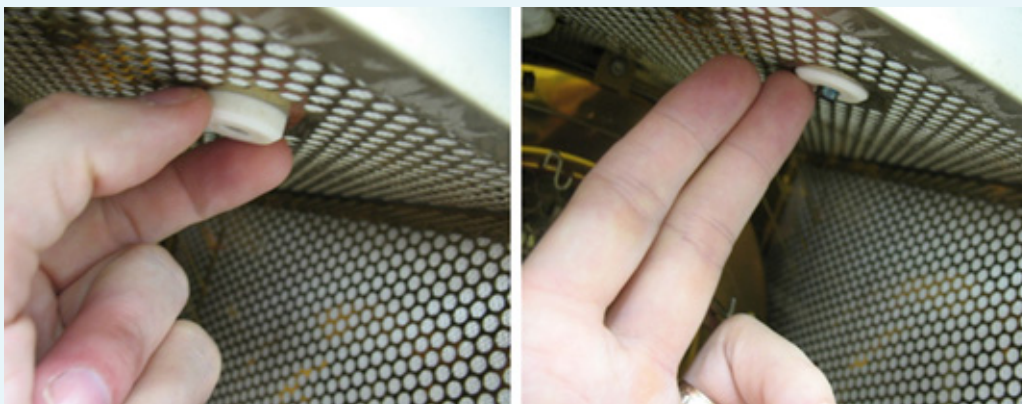

This will expose the end of the insert and make it easier to locate the column in the central hole.

After the final trimming and installation of the column, make sure to pull the insulator back to its lowest position.

## Checking for Leaks

A hand held helium leak detector may be used to check for helium leaks around the inlet and column connection. With the helium conservation enabled, approximately 4 mL/min of helium will enter the injection port.


For maximum sensitivity of the leak detector, set the split flow to **Off** during leak check. This results only in column and septum purge flows.



The composition of the gas mixture in the inlet will be mostly helium and allow good sensitivity. For bulk nitrogen leak detection, the GC column exit can be plugged, and the automated pressure drop leak detection of the GC used for determining the pressure drop.

The helium delivery block should be parked in the rear “blocked” position of [Figure 475](#) if leak checking is done based on pressure drop. The block must be replaced in the front position to resume operation.


**Figure 475.** Removing the Helium Tube from the Gas Delivery Block

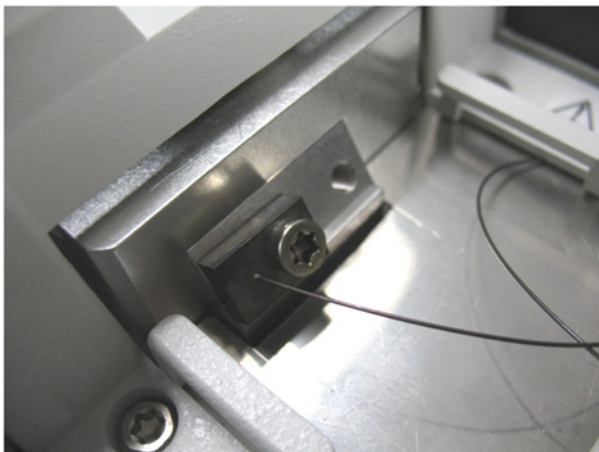

## Adding a PDD Module

This section provides instructions for adding a front/back PDD detector module. According to the configuration of your TRACE 1300/TRACE 1310, a dummy module is present into the free site where the detector module is not installed. See [Figure 476](#).

**Figure 476.** Add a Front/Back Detector Module

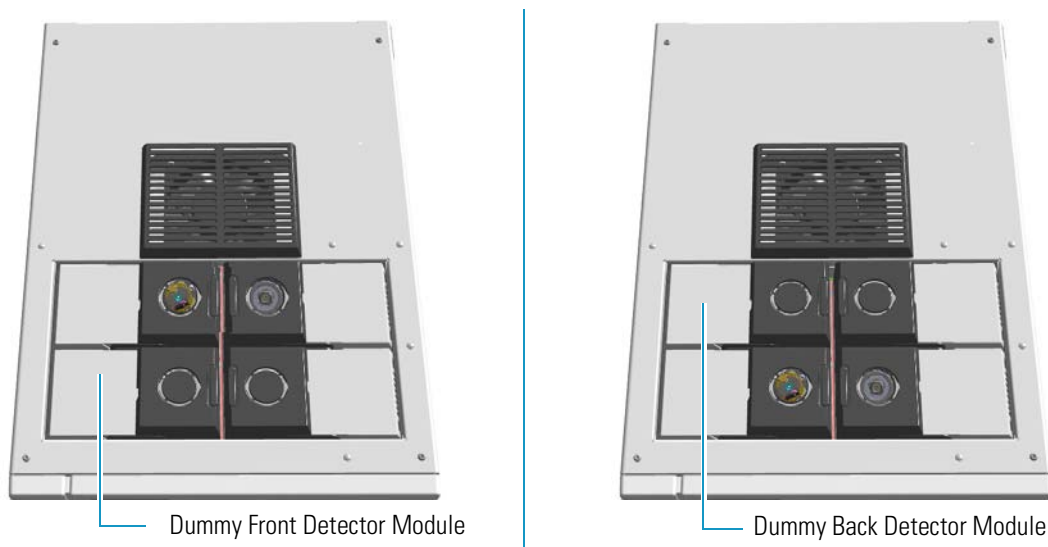

**ATTENTION** Where a dummy module is installed, the gas connection is blocked by a plug.

## Getting Started

To properly installing and connecting the PDD module see sequentially the following sections:

- “Installing a PDD Module” on [page 459](#)
- “Plumbing the Gas Lines” on [page 463](#)
  - “Installing and Purging the Gas Regulator” on [page 464](#)
  - “Installing and Purging the Helium Purifier” on [page 464](#)
  - “Connecting the Discharge Gas Supply” on [page 465](#)
- “Connecting the High Voltage and Pulses Cables” on [page 465](#)
- “Installing the Capillary Column” on [page 470](#)
- “Installing the Packed Column” on [page 470](#)
- “Testing for Leaks” on [page 470](#)
- “Performing Initial Power Up” on [page 471](#)




## Installing a PDD Module

### ❖ To install a PDD detector module

The result of the operation is shown in [Figure 477](#).

**Figure 477.** Result of the PDD Installation

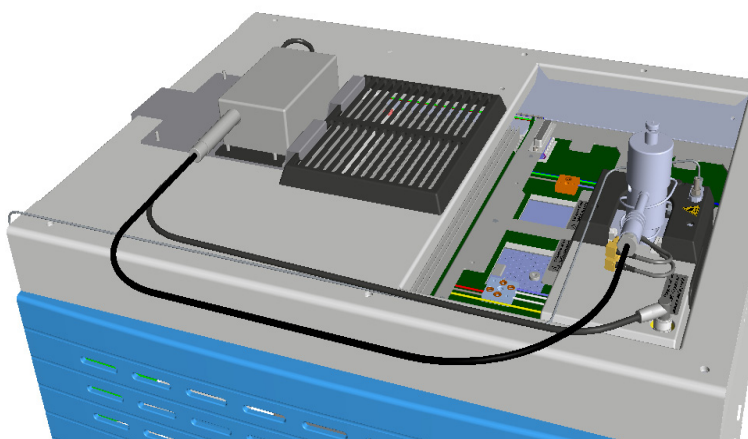

**Note** The installation of a single PDD detector is considered. The injector modules are deliberately missing for convenience.

1. Put the GC in standby condition.
2. Cool the oven, injector, and detector to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
6. Remove the dummy module from the position where the detector module will be installed.
  - a. Open the module flap cover.
  - b. Using a T20 Torxhead screwdriver, unscrew and remove the two captive fixing screws.

- c. Keeping the dummy module flap cover open, lift up the module from its seat in the injector/detector housing. Place the dummy module on a clean surface.
- d. Keep the gas block plug of the gas connections installed. See [Figure 447](#).


**Figure 478.** Detector Gas Block Plug



7. Open the front door of the GC.
8. Plug and fix the detector module into the main frame.
  - a. Place the module in its seat. Be sure to insert the 25-pin male connector, on the bottom of the module, into the 25-pin female connector on the detector seat of the detector housing.

**Note** Because the module flap cover of the PDD detector cannot be opened, two holes are provided on the module flap cover for accessing the captive fixing screws. See [Figure 479](#).

**Figure 479.** Holes on the Module Flap Cover

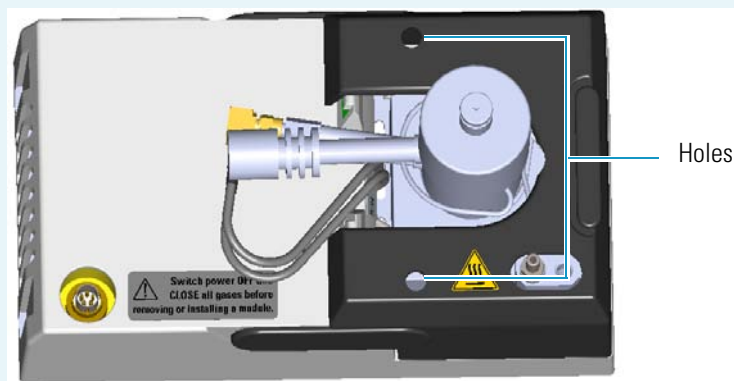

- b. Insert a T20 Torxhead screwdriver into the holes provided on the module flap cover and tighten the captive fixing screws without overtightening.




**ATTENTION** To maintain the correct alignment the screws must be tightened in turn. Tighten each screw only a small amount before moving to the next screw. Repeat until all are secure.

9. Mount the support bracket of the pulse generator. See [Figure 480](#).


**Note** The support bracket can hold two pulse generators for supplying two PDD detectors on the GC.

**Figure 480.** Support Bracket

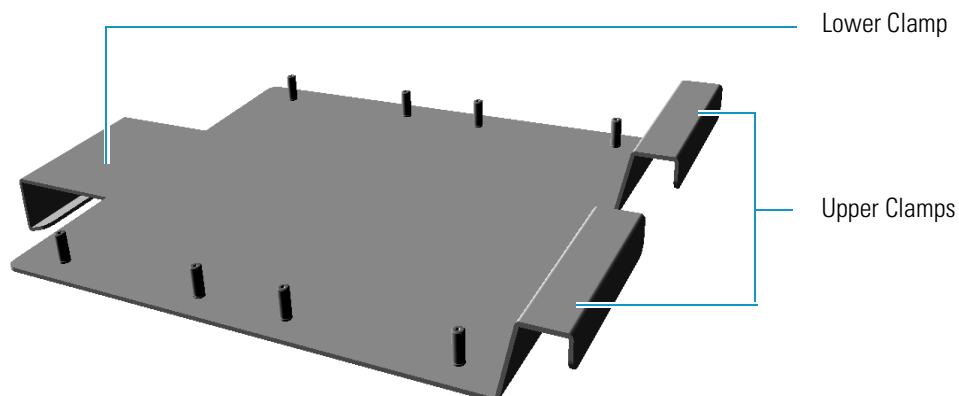


- a. Insert the upper clamps of the support bracket into the last slots of the cooling fan grid. See [Figure 481](#).

**Figure 481.** Support Bracket Installation (1)

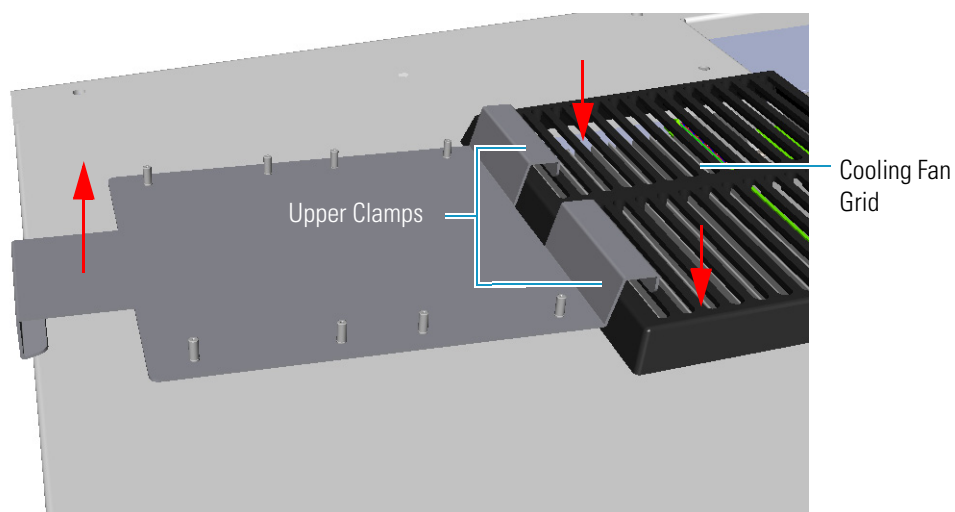

- b. Lift the bottom part of the support bracket and push it back until the lower clamp hooks the rear of the GC top cover. See [Figure 482](#).

**Figure 482.** Support Bracket Installation (1)

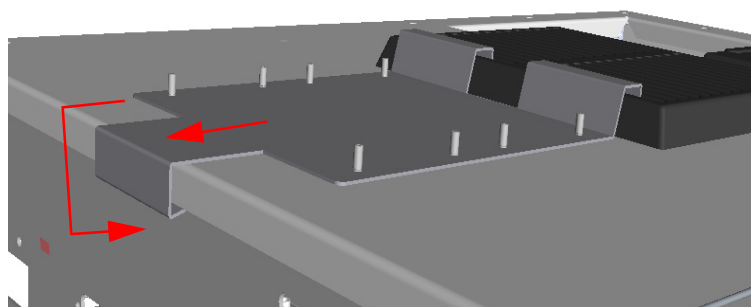

**Note** Two support brackets can be installed for supplying a maximum of four PDD detectors, two installed on the GC and two installed on the TRACE 1310 Auxiliary Oven. See the example of [Figure 483](#) and [Figure 484](#).

**Figure 483.** Installation of two Support Bracket




**Figure 484.** Installation of Four PDD Detectors



10. Place the pulse generator on the support bracket.
  - a. Place the pulse generator on a convenient position of the support bracket by aligning the four holes on the base of the generator with the corresponding four rivets on the support bracket as shown in the example of [Figure 485](#).

**Figure 485.** Pulse Generator Installation

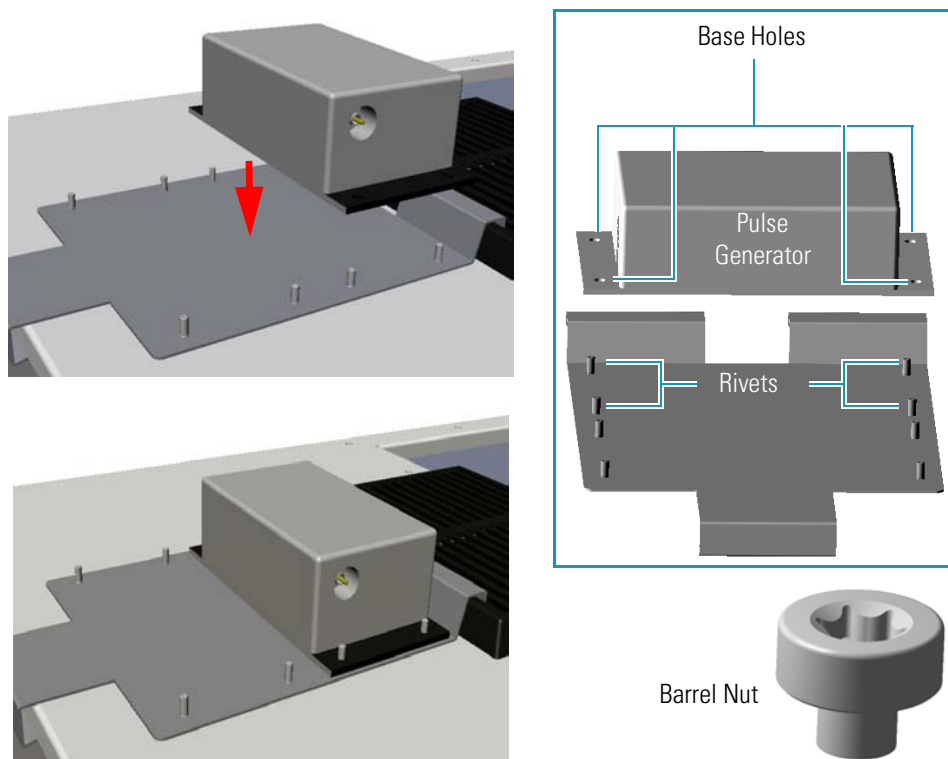


- b. Fix the pulse generator using the four barrel nuts provided. Tighten the nuts using a T30 Torxhead screwdriver.

## Plumbing the Gas Lines

Figure 486 shows the diagram of the gas lines connections of the detector system.

**Figure 486.** Gas Line Connections

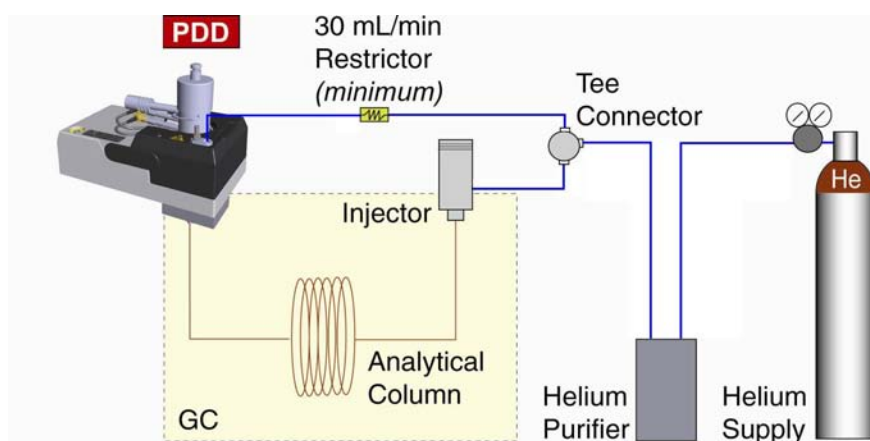

Before connecting the gas lines verify that:

- The pressure regulators are commercial ultra-pure grade regulators with stainless steel diaphragms. Regulators with diaphragms made of neoprene or others elastomers should never be used.
- The connecting tubes are thoroughly cleaned and baked before use.
- The gas regulator and the helium purifier must be properly purged. See the following sections for further details: “[Installing and Purging the Gas Regulator](#)” on page 464, “[Installing and Purging the Helium Purifier](#)” on page 464, and “[Connecting the Discharge Gas Supply](#)” on page 465.



**WARNING** Never use leak detecting fluids on any part of the system.

## Installing and Purging the Gas Regulator



**WARNING** To avoid injury, before starting make sure to strictly respect all the Regulations concerning the manual handling of gas cylinders under pressure.

### ❖ To install and purge the gas regulator

1. Install the pressure regulator on the cylinder.
2. Making sure the regulator is closed, open the cylinder allowing the regulator to pressurize.
3. Close the cylinder and check the regulator for pressure leaks. Observe the needle of the regulator output pressure gauge for 15 minutes.
4. Open the regulator allowing the pressure to relief.
5. Repeat five times from step 2.

## Installing and Purging the Helium Purifier

### ❖ To install and purge the helium purifier

Before starting, make sure the helium purifier outlet tube is capped.

1. Connect the helium cylinder pressure regulator to the inlet port of the helium purifier by using the appropriate connecting tube and fittings.
2. Turn the output pressure regulating knob clockwise until the gauge registers 345 kPa (50 psig).
3. Wait five minutes for equilibrium, then turn the regulating knob all the way counterclockwise.


4. Observe the needle of the output pressure gauge for 15 minutes. There will be a slight initial drop. If it does not move after that, consider all the connections are tight.
5. If necessary, use an electronic leak detector to locate any leaks. If a leak detector is not available, tighten all the fitting and repressurize the system for another test.
6. Uncap the outlet tube of the helium purifier and purge the system for 15 to 30 minutes at 60 - 80 mL/min to eliminate air from the purifier getting material.

## Connecting the Discharge Gas Supply

### ❖ To connect the discharge gas supply

1. Connect the helium purifier outlet port to a port of the Tee connector provided by using the 1/16-in. o.d. connecting tube provided.
2. Connect the second port of the Tee connector to the carrier gas inlet port, located on the rear panel of the GC, by using a sufficient piece of the stainless steel connecting tube provided and the appropriate fitting.
3. Connect the last port of the Tee connector to the inlet of the discharge gas restrictor.
4. Connect the outlet tube of the restrictor to a flow measuring device and adjust the helium pressure to obtain a flow of about 30 mL/min.
5. After setting the flow rate, connect the outlet of the restrictor to the discharge gas inlet on the PDD module. See [Figure 487](#).

**Figure 487.** Gas Supply Connections

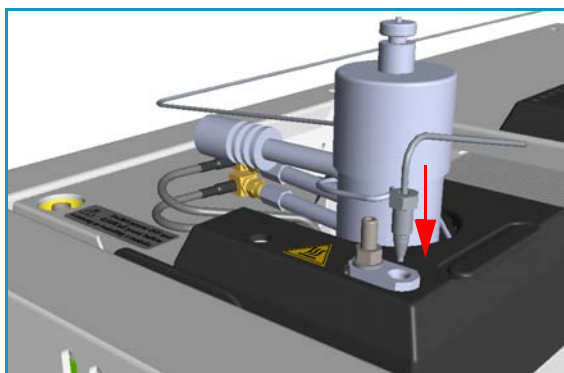

## Connecting the High Voltage and Pulses Cables


### ❖ To connect the cables from the pulse generator to the PDD module



**WARNING** Never the high voltage discharge cable must be disconnected when the pulse generator is turned on. Dangerous high voltage is present: 3700 V Peak; 170 mA Peak.


To avoid accidental disconnection of the high voltage discharge cable (mainly) and the pulses cable from the pulse generator, the safety covers must be installed. See [Figure 488](#).

**Figure 488.** Pulse Generator Safety Covers

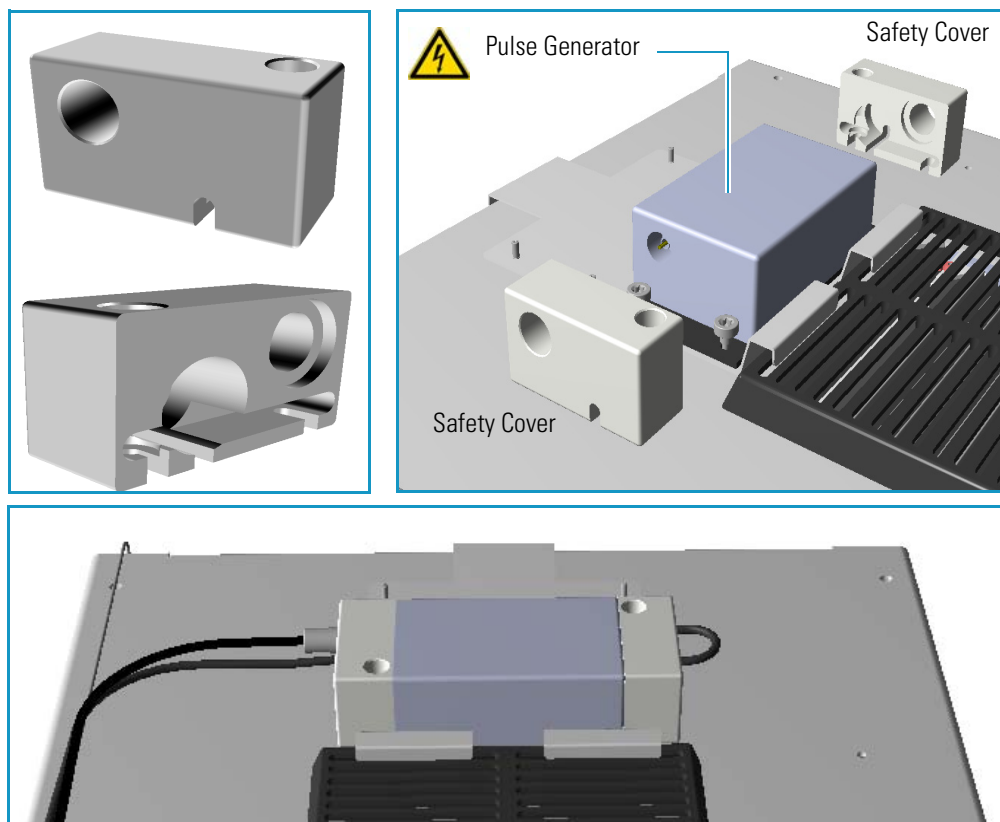

1. Connect the high voltage discharge cable.
  - a. On the high voltage side, loosen the right nut that fix the pulse generator to the bracket. Do not loosen the left nut, keep it tightened. Place the safety cover until the nuts couple the slots on the bottom of the safety cover. See [Figure 489](#).


**Figure 489.** High Voltage Discharge Cable Connection (1)

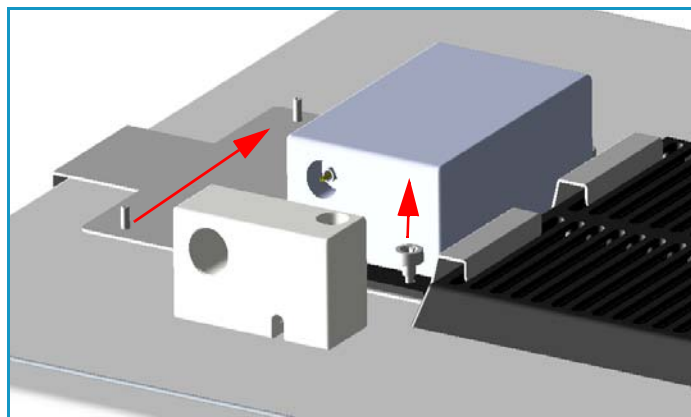

- b. Tighten the right nut using a T30 torxhead screwdriver. See [Figure 490](#).




**Figure 490.** High Voltage Discharge Cable Connection (2)

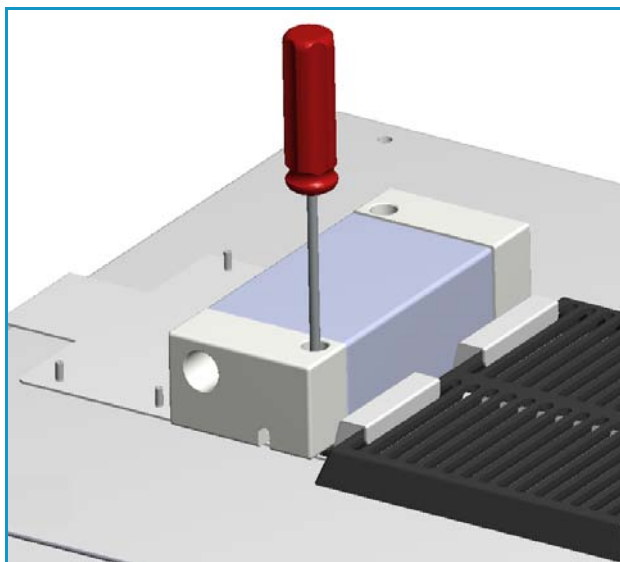

- c. Insert the connector of the high voltage discharge cable, coming from the high voltage discharge electrode on the PDD module, into the connection provided in the pulse generator passing through the hole on the safety cover. See [Figure 491](#).


**Figure 491.** High Voltage Discharge Cable Connection (3)

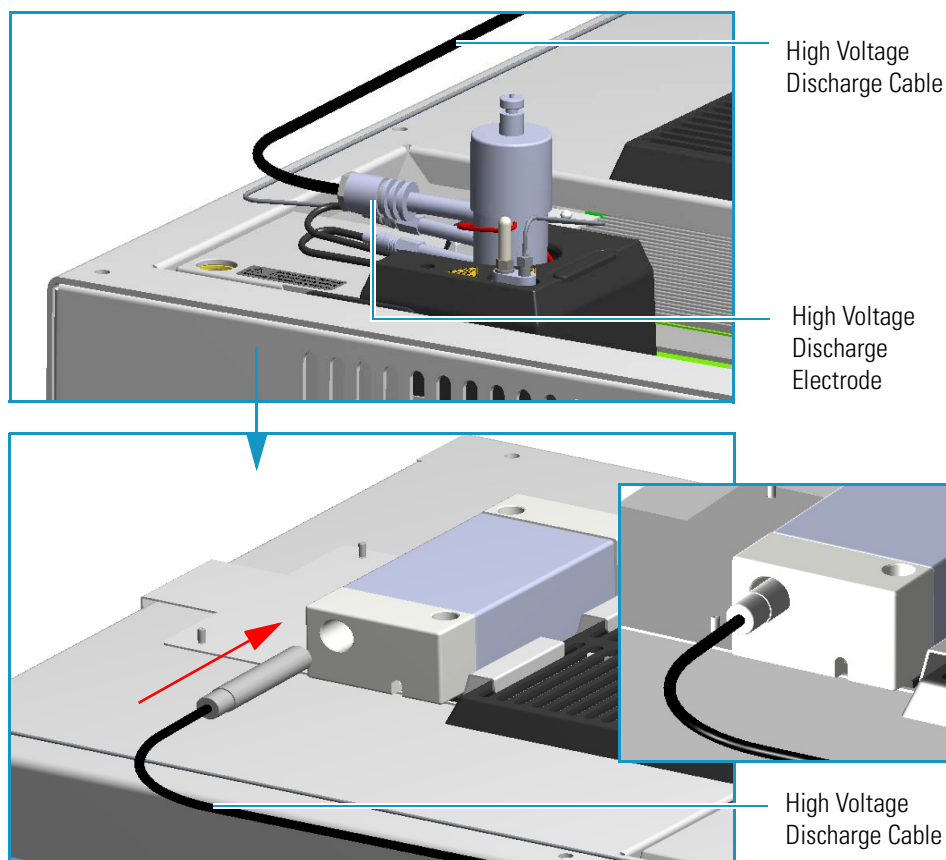

2. Connect the pulses cable.


- a. Plug the 4-pin pulses cable connect into the 4-pin connector on the pulse generator. See [Figure 492](#).

**Figure 492.** Pulses Cable Connection (1)

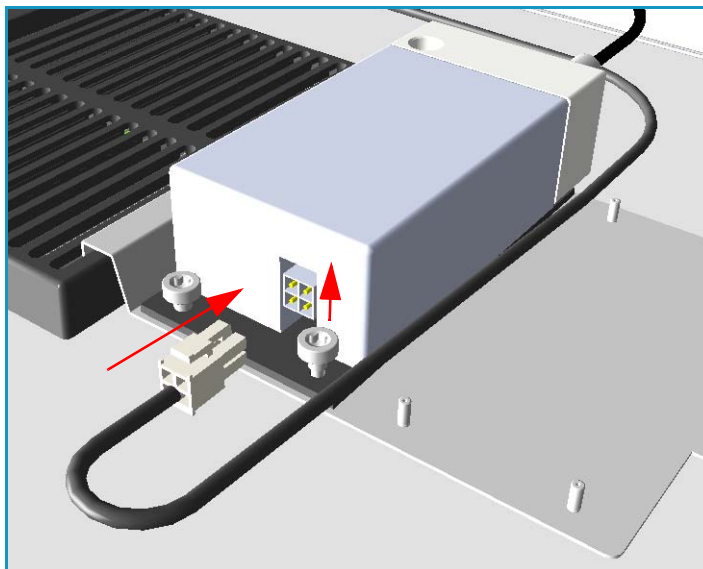


- b. On the pulse generator control side, loosen the right nut that fix the pulse generator to the bracket. Do not loosen the left nut, keep it tightened. Place the safety cover until the nuts couple the slots on the bottom of the safety cover paying attention to guide the cable into the slot located on the bottom of the safety cover. See [Figure 493](#).

**Figure 493.** Pulses Cable Connection (2)

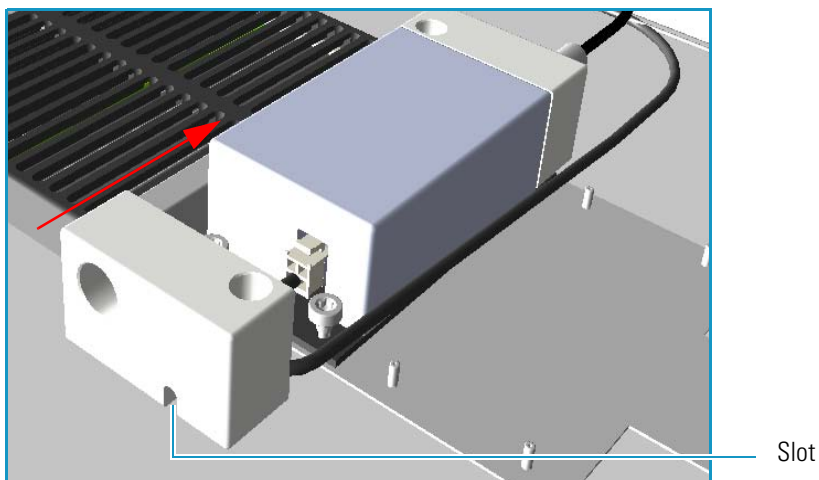


- c. Tighten the right nut using a T30 torxhead screwdriver. See [Figure 494](#).

**Figure 494.** Pulses Cable Connection (3)

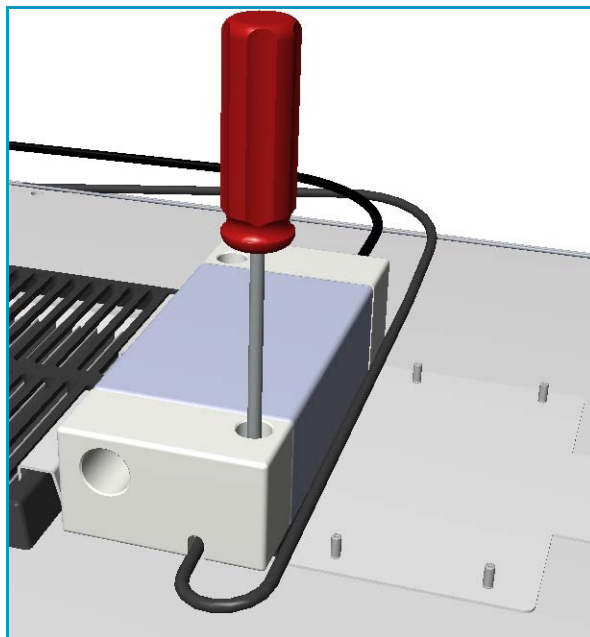

- d. Connect the cable between the pulse generator and the pulse generator control on the PDD module. See [Figure 495](#).

**Figure 495.** Pulses Cable Connection (4)

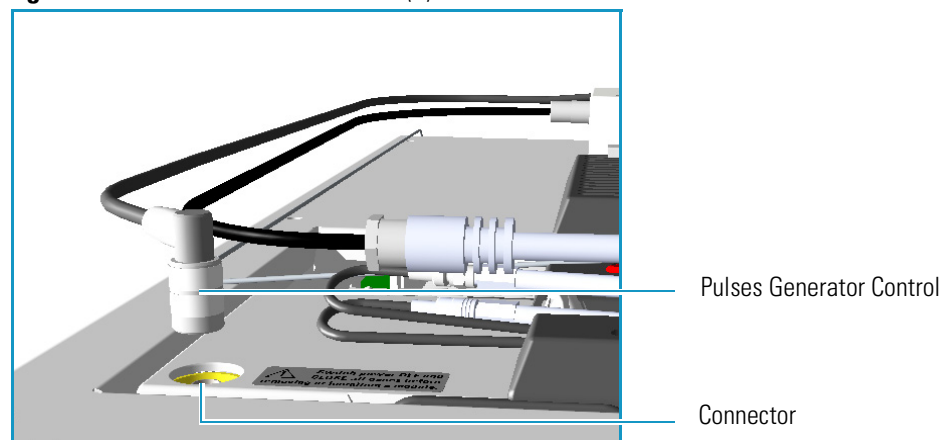

**WARNING** If for any reasons it is necessary to disconnect the cables from the pulse generator, the safety covers must be removed proceeding in the reverse order in which they was installed. Before starting, make sure that the GC is powered Off.

## Installing the Capillary Column

For PDD, the column must penetrate **136 mm** inside the pre-installed capillary column adapter.

### ❖ To install the capillary column

1. Make a mark on the column 136 mm from the end.
2. Remove the knurled nut column inlet at the bottom of the detector. Slide the nut over the end of the column, followed by the appropriate column ferrule.
3. Seat the ferrule in the detail of the column adapter and begin sliding the column through the capillary column adapter and into the column inlet.
4. Get the nut started on the threads and tighten it until you feel it contact the ferrule, then back off half a turn.
5. Slide the column into the column inlet until the mark is flush with the surface of the knurled nut, and secure the column in the adapter by tightening the knurled nut finger tight only.

## Installing the Packed Column

### ❖ To install the packed column

1. Remove the pre-installed capillary column adapter.
2. Replace pre-installed capillary column adapter with the **packed columns adapter** that penetrates into the PDD cell for the correct length.
3. Connect the packed column to the packed column adapter.

## Testing for Leaks

It is critical for the system to be leak-tight. Leak test is strongly recommended before operating with PDD.

### ❖ To test for leaks

1. Open the discharge gas supply.
2. Cap the tube and pressurize the entire system with helium to 138 kPa (20 psi).
3. If the system does not hold pressure, check all the fittings with an electronic helium leak detector. DO NOT use leak detecting liquids.
4. Tighten fittings as required.

## Performing Initial Power Up

### ❖ To perform the initial power up



**WARNING** During normal operation, the detector produce ultraviolet energy (UVA, UVB), some of which may be emitted. Do not watch the arc without eye protection.

1. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
2. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker) to the position I.
3. Close the front door of the GC.
4. Set the detector temperature to 100 °C and allow time for the detector to reach the set temperature, and for the helium purifier to reach the optimum working condition.
5. Turn on the pulse generator.
6. Check the standing/background current. Optimum range is 600 - 2000 pA at 100°C. Lower current indicates a clean, leak-free system.
7. The recommended detector temperature is 20 °C above the column temperature, with a minimum of 100 °C. Set the detector to the operating temperature required for the intended analysis. When the detector has reached the set temperature, read and record the standing current.
8. Start carrier flow, then read the standing current. The difference between this reading and the one previous is the ionization of the combined impurities in and eluting with the carrier gas. The smaller the difference, the better the quality of the gas exiting the column.
9. Set the column oven to the temperature required for the intended analysis. When the oven reaches the set temperature, read the standing current. The difference with the previous reading is the ionization of the column bleed. The smaller the difference, the better the column is conditioned.

From this point, the standing current should be observed and logged after any system change. In addition, logging the standing current (with and without the column) on a regular basis is an effective monitor of system integrity (leak tightness and cleanliness).

We also recommend tracking the internal standard (quantity on column/area count) for sensitivity continuity.

## Adding a Generic Detector Interface

This section provides instructions for installing the detector and control modules of the Generic Detector Interface on your TRACE 1300/1310 GC or TRACE 1310 Auxiliary Oven.



**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.

### Preliminary Operations

Before starting, the following preliminary operation must be carried out.

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.



**CAUTION** If the GC is coupled to a TRACE 1310 Auxiliary Oven, unplug the power cable from the AC Input connector into the back of the TRACE 1310 Auxiliary Oven and from the wall outlet.

4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.

### Getting Started

For properly installing the Generic Detector Interface see the following sections:

- “Removing the GC Back Cover” on [page 473](#)
- “Assembling the GDI Electrical Interface” on [page 474](#)
- “Replacing the Encapsulated Flow Restrictors” on [page 476](#)
- “Installing and Connecting the GDI Electrical Interface” on [page 477](#)
- “Connecting a GDI Interface to the TRACE 1310 Auxiliary Oven” on [page 481](#)

- “Installing a GDI Mechanical Module” on page 482
- “Connecting the Detector Gas Tubings to the Manifolds” on page 484
- “Connecting Heater and Signal Cables” on page 488
- “Restarting the GC” on page 489
- “Performing the Third-party Detector Start-up and Optimization” on page 489
- “Configuring and Setting GDI Detector” on page 490


## Removing the GC Back Cover

Remove the back cover for accessing the compartment on the back of the GC.

### ❖ To remove the back cover

1. Remove the cover.
  - a. Using a T20 Torxhead screwdriver, unscrew the four screws that secure the back cover to the GC. See Figure 496.

**Figure 496.** Back Cover Removal

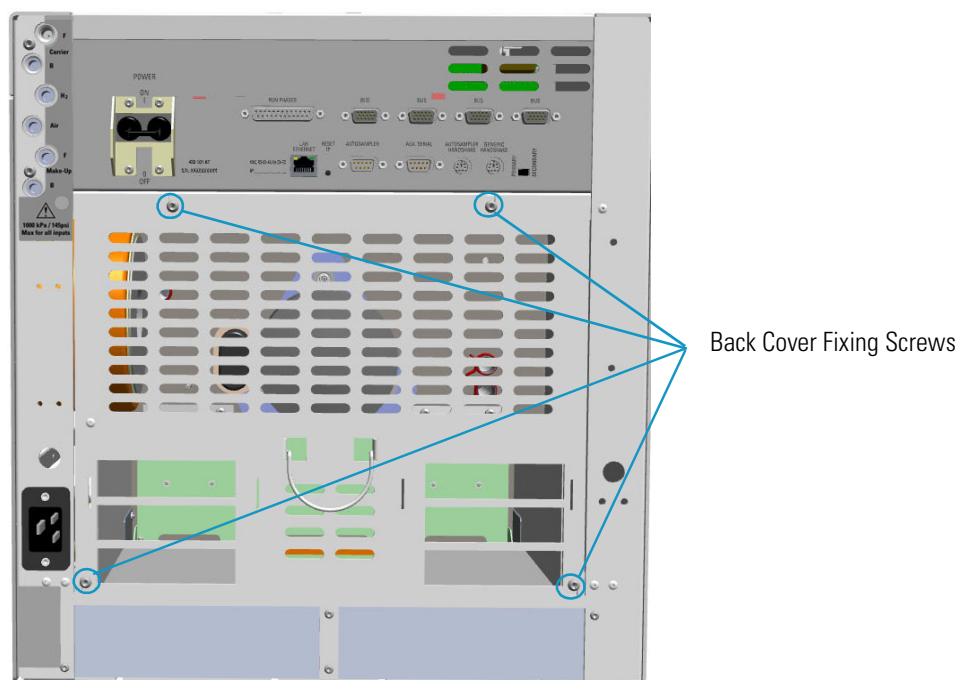


- b. Lift the cover off using the cover handle. Be aware that the ground wire is attached to the back cover.

**Note** Pay attention to the positioning of the ground wire plug, so it can be reconnected in the same way it was removed.

## Assembling the GDI Electrical Interface

If not already done, you must connect the detector cable and its cable grommet into the GDI electrical interface. See [Figure 497](#).


**Figure 497.** GDI Electrical Interface Back View

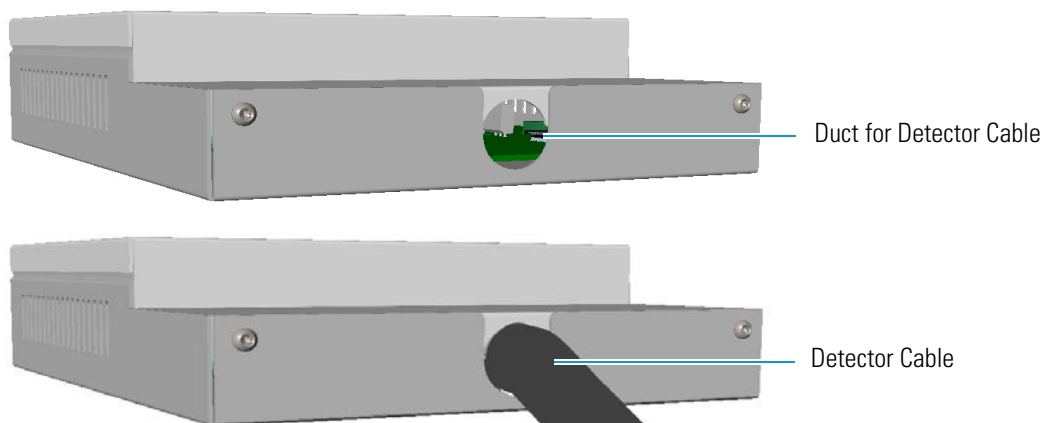

### ❖ To assemble the GDI electrical interface

1. Unscrew the fixing screws of the GDI electrical interface using a T10 torxhead screwdriver, and remove the top cover. See [Figure 498](#).


**Figure 498.** Assemble GDI Electrical Interface (1)

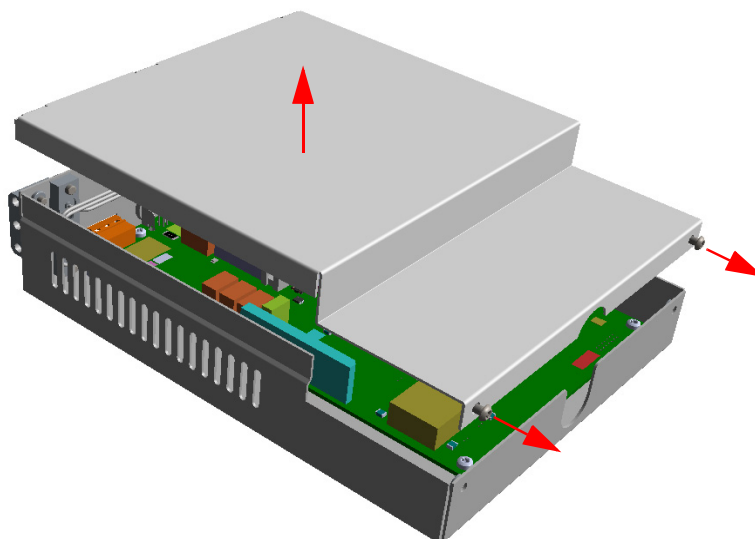

2. On the board inside the control module look for the 26-pin male connector marked **J4**. See [Figure 499](#).




**Figure 499.** Assemble GDI Electrical Interface (2)

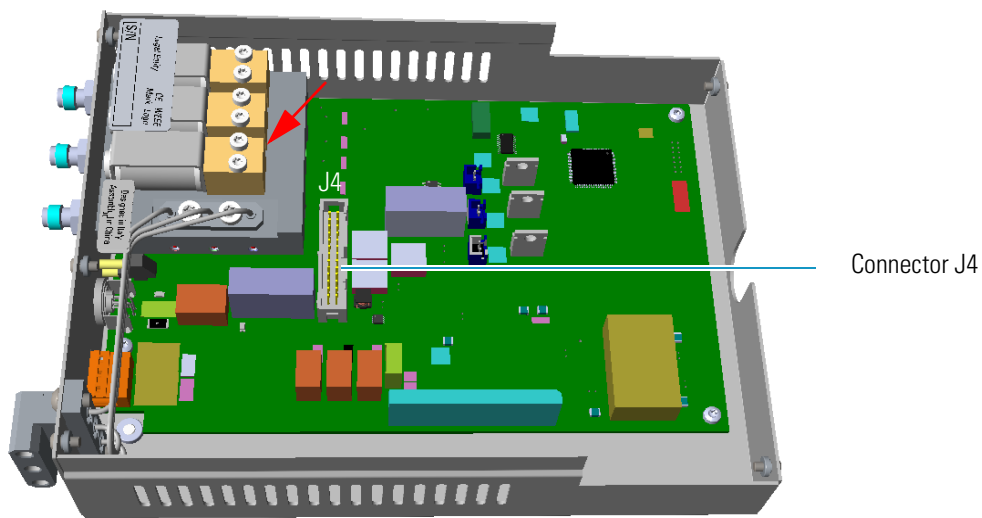

3. Connect the 26-pin female connector of the detector cable to the 26-pin male connector marked **J4** on the board.
4. Guide the detector cable and its cable grommet out of the control module through the duct. See [Figure 500](#).

**Figure 500.** Assemble GDI Electrical Interface (3)

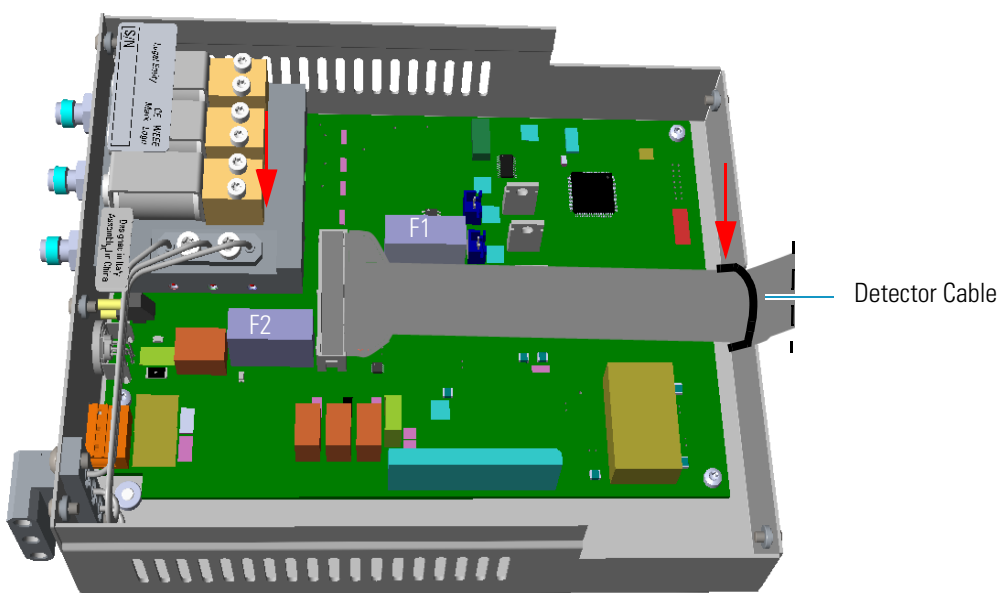

**ATTENTION** Two protecting fuses **F1** and **F2** are present on the electronic board of the GDI electrical interface for the protection of the **24 Vdc** and **48 Vac** respectively:


- **F1** = T300mA 250V; (5 x 20 mm)
- **F2** = T2A 250V; (5 x 20 mm)

5. Remount and fix the top cover by using the fixing screws.

## Replacing the Encapsulated Flow Restrictors

This section provides instruction for replacing the 1/4-in. OD x 1/4-in. length encapsulated flow restrictors. The flow restrictors are located inside the manifold marked **Gas Outputs** on the GDI electrical interface. See [Figure 501](#).


**Figure 501.** Gas Outputs Manifold

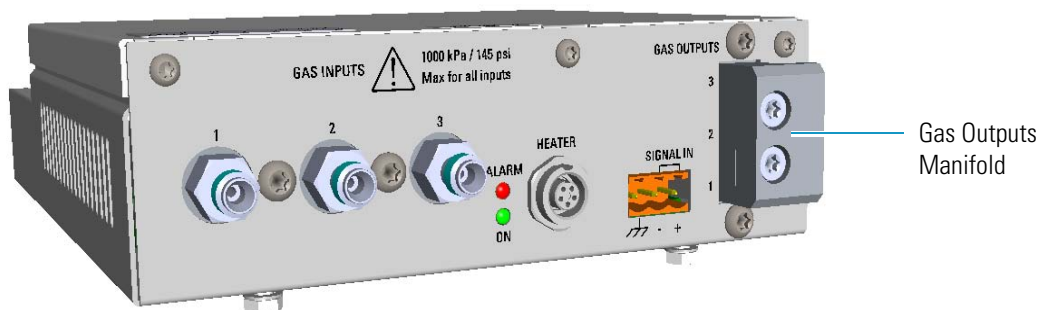

### ❖ To replace an encapsulated flow restrictor


1. Remove the Gas Outputs manifold from the GDI electrical interface by unscrewing its fixing screws using a T20 Torxhead screwdriver. See [Figure 502](#).

**Figure 502.** Gas Outputs Manifold

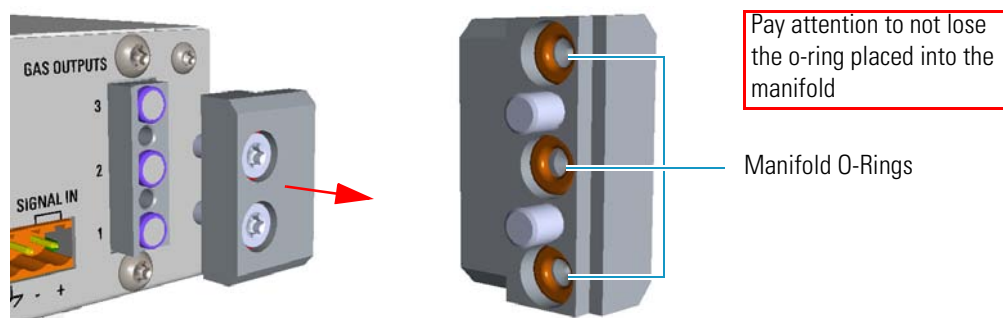

2. Remove the encapsulated flow restrictor of interest from its seat using tweezers. See [Figure 503](#).

**Figure 503.** Flow Restrictor Replacement

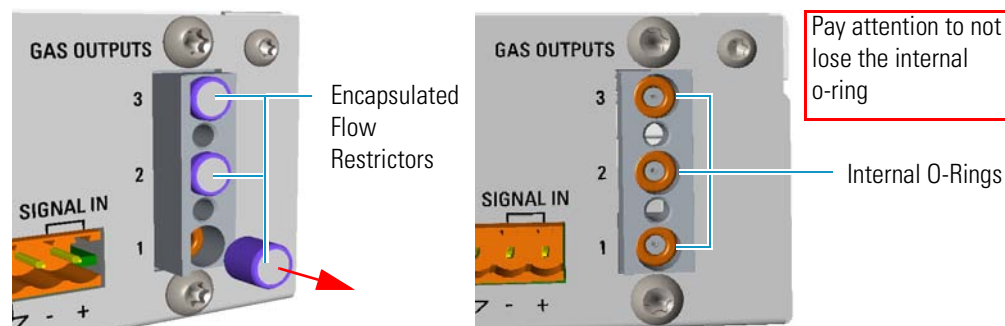

3. Replace the encapsulated flow restrictor with the ones required by the third-party detector, then reinstall the manifold proceeding in the reverse order in which it was removed.
4. Calibrate the new full scale of the flow restrictor.
  - a. In the **GDI Configuration** page, set gas type to the real gas in use and full scale to 100 mL/min whatever restrictor is installed.  
  
See the section “[Configuring and Setting GDI Detector](#)” on [page 490](#), and refer to Chapter 2 and Chapter 4 of the *TRACE 1300/TRACE 1310 User Guide*.
  - b. Set relevant channel flow rate in the **Method/Instrument Control** page to the full scale value (100 mL/min). Make sure the gas input pressure is at least 60 psig (414 kPa), and verify that the actual flow displayed by the GC reaches 100 mL/min as set.  
  
See the section “[Configuring and Setting GDI Detector](#)” on [page 490](#), and refer to Chapter 2 and Chapter 4 of the *TRACE 1300/TRACE 1310 User Guide*.
  - c. Check the real flow rate with an external flowmeter.
  - d. Use the real reading of the flowmeter as the new full scale, and set it in the **GDI Configuration** page.

## Installing and Connecting the GDI Electrical Interface

This section provides the instruction for installing and connecting the GDI electrical interface.



**CAUTION** It is mandatory that the GDI electrical interface must be placed into an housing on the back of the GC. If both the housings are already occupied by other external modules, one of these must be removed and placed beside the GC, and in the most comfortable position for the user to have free access to the connections.

### ❖ To install and connect the GDI electrical interface

1. Remove the cover of the external modules housing where installing the module.  
See [Figure 504](#).

**Figure 504.** Housing Cover Removal




## 9 Adding Modules

### Adding a Generic Detector Interface


- a. Using a T20 Torxhead screwdriver, unscrew and remove the left or right housing cover screws, then remove the cover from the housing.
2. Install the GDI electrical interface into the housing
  - a. Loosen the two hexagonal screws under the module. See [Figure 505](#).

**Figure 505.** GDI Module Installation (1)

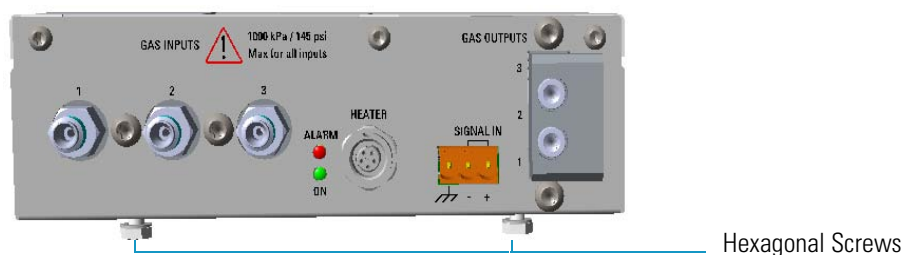

- b. Carefully place the GDI electrical interface into the left or right housing. Guide the detector cable into the electronic compartment of the GC, next push the module until the hexagonal screws couple with the slots on the floor of the GC. See [Figure 506](#).

**Figure 506.** GDI Module Installation (2)

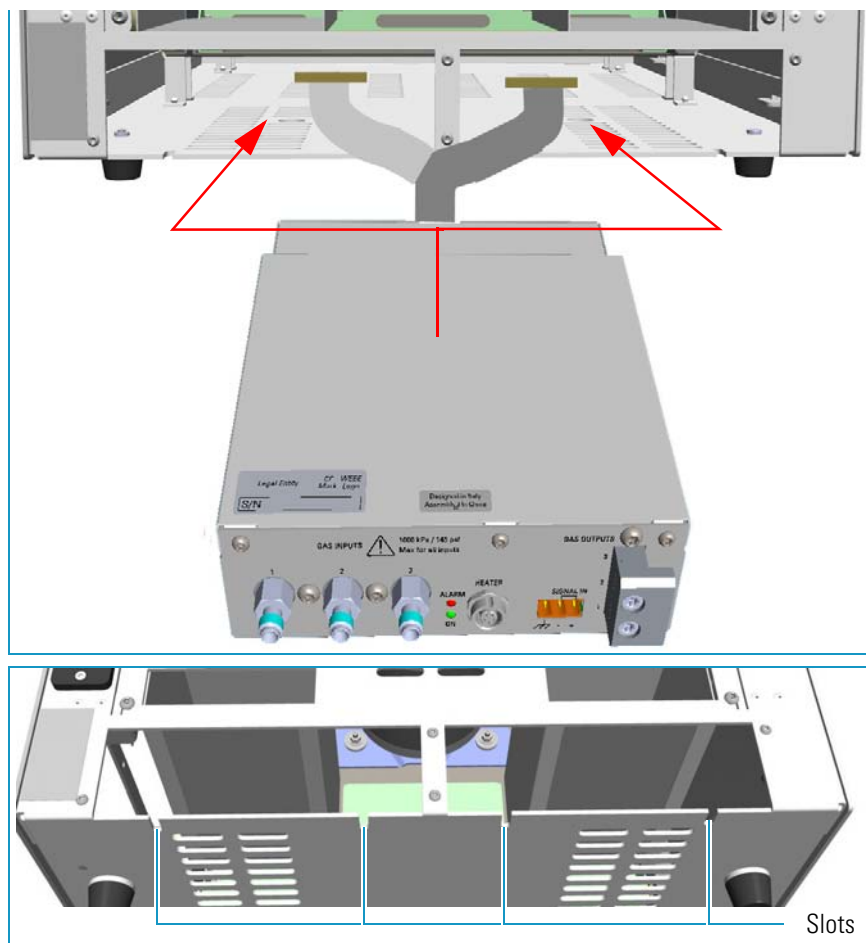

- c. Finger-tighten the hexagonal screws slightly, or use a 10-mm wrench.


**Note** Always keep the hexagonal screws in their place. This allows you an easier removal of the GDI electrical interface when necessary.



**ATTENTION** If an GDI mechanical module will be installed on the TRACE 1300/TRACE 1310 GC, continue the procedure from the [step 3](#). If an GDI mechanical module will be installed on the TRACE 1310 Auxiliary Oven, jump to the section “Connecting a GDI Interface to the TRACE 1310 Auxiliary Oven” on [page 481](#).

3. Connect the 26-pin connector of the detector cable coming from the GDI electrical interface to the backplane board in the electronic compartment.
  - a. On the backplane board disconnect the flat cable from the connector marked **J13 DET.FRONT** if the GDI mechanical module is installed in the **Front** site, or from the connector marked **J5 DET.REAR** if the GDI mechanical module is installed in the **Back** site. See [Figure 507](#).

**Figure 507.** GC Backplane Board Layout

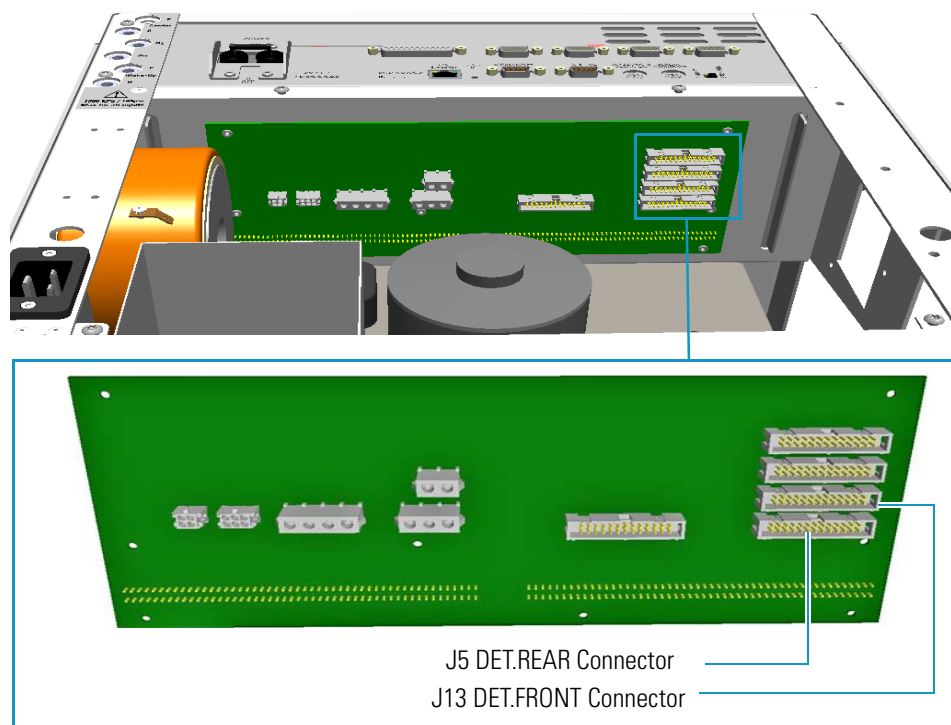

- b. Connect the detector cable to the connector **J13 DET.FRONT**, or **J5 DET.BACK** accordingly.

#### **ATTENTION** TRACE 13100/1310 coupled with a TRACE 1310 Auxiliary Oven

If you are installing the GDI mechanical module in the **Front** or **Back** site of a TRACE 1300/1310 GC coupled with a TRACE 1310 Auxiliary Oven, the detector cable of the GDI electrical interface must be connected as follows:




- **Mechanical module installed in the Front site:** On the backplane board disconnect the flat cable from the connector marked **J13 DET.FRONT**. Connect the detector cable to the connector **J13 DET.FRONT**.
- **Mechanical module installed in the Back site:** Do NOT DISCONNECT the flat cable from connector marked **J5 DET.REAR**. Connect the detector cable directly to the additional connector on the flat cable.

4. Reconnect the ground wire to the back cover terminal.
5. Replace the back cover of the GC proceeding in the reverse order in which it was removed.

The result of the installation is shown in [Figure 508](#).

**Figure 508.** GDI Electrical Interface Installed into the GC

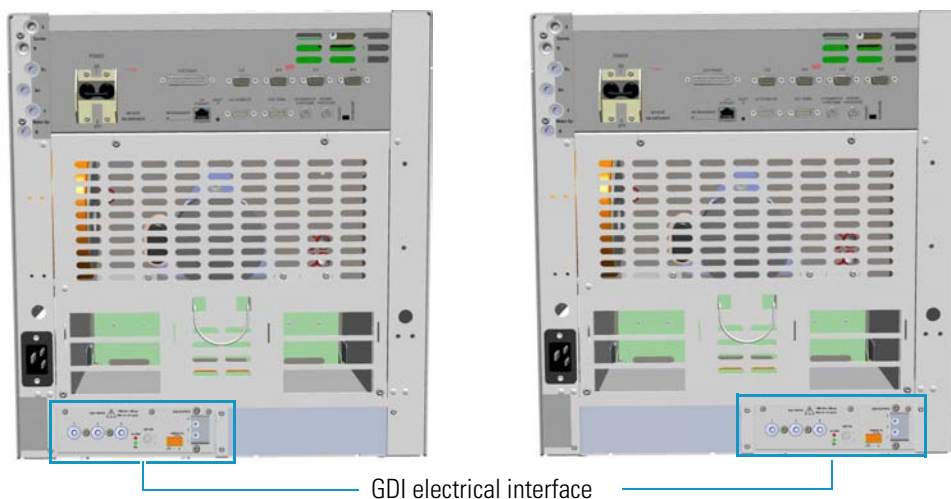

**Note** Continue the procedure of installation following the instruction reported in the section “[Installing a GDI Mechanical Module](#)” on [page 482](#).

## Connecting a GDI Interface to the TRACE 1310 Auxiliary Oven

The GDI mechanical module can be installed as auxiliary detector into the position **Aux L** or **Aux R** of the auxiliary detector assembly.

The detector cable of the GDI electrical interface, installed into an external module housing on the back of the GC, must be connected to the VOBP-HRM board located into the electronic compartment of the TRACE 1310 Auxiliary Oven. See the following procedure.

### ❖ To connect the GDI electrical interface to the TRACE 1310 auxiliary oven



**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.

Before starting make sure that the TRACE 1300/1310 GC and the TRACE 1310 Auxiliary Oven are powered off and the power cables disconnected.

1. Carefully decouple the TRACE 1310 Auxiliary Oven from the GC for creating enough operating space for guiding the detector cable from the GDI electrical interface through the GC into the electronic compartment of the TRACE 1310 Auxiliary Oven.



**WARNING** Before proceeding this operation make sure to disconnect properly the inner tube. Refer to the section Coupling the TRACE 1310 Auxiliary Oven to the GC on the TRACE 1310 Auxiliary Oven Instruction Manual, and operating in the reverse order in which it was mounted.

2. In the electronic compartment look for the VOBP-HRM board. It is located over the power section of the TRACE 1310 Auxiliary Oven. See [Figure 509](#).


**Figure 509.** TRACE 1310 Auxiliary Oven: VOBP-HRM Board



3. Connect the 26-pin connector of the detector cable coming from the GDI electrical interface to the VOBP-HRM board in the electronic compartment.
  - a. On the board VOBP-HRM disconnect the flat cable from the connector marked **J29 AUX DET LEFT** if the GDI mechanical module is installed in the **Left** site, or from the connector marked **J13 AUX DET RIGHT** if the GDI mechanical module is installed in the **Right** site. See [Figure 510](#).



**Figure 510.** VOBO-HRM Board Layout

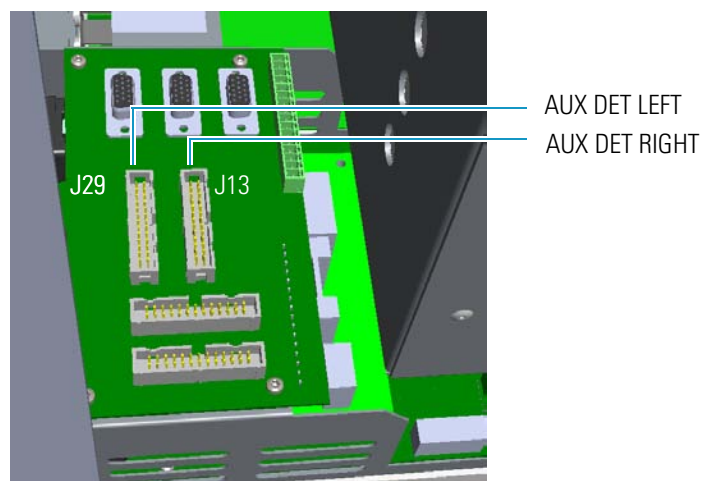


- b. Connect the detector cable to the connector **J29 AUX DET LEFT**, or **J13 AUX DET RIGHT** accordingly.
4. Re-couple the TRACE 1310 Auxiliary Oven to the GC.

**Note** Continue the procedure of installation following the instruction reported in the section “[Installing a GDI Mechanical Module](#)” on [page 482](#).

## Installing a GDI Mechanical Module

This section provides instructions for adding a **Front/Back/Aux L/Aux R** GDI mechanical module. According to the configuration of your TRACE 1300/TRACE 1310 GC, a dummy module is present into the free site where a detector module is not installed. See [Figure 446](#).

**Figure 511.** Add a Front/Back/Aux L/Aux R Detector Module

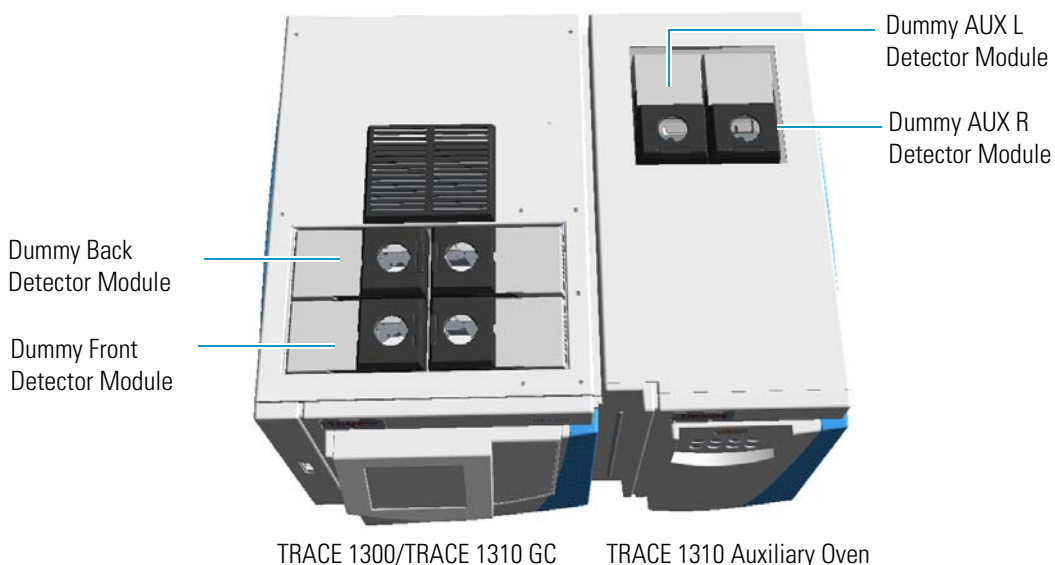






**ATTENTION** Where a dummy module is installed, the gas connection is blocked by a plug.

❖ **To install a GDI mechanical module**


1. Remove the dummy module from the position where the detector module will be installed.
  - a. Open the module flap cover.
  - b. Using a T20 Torxhead screwdriver, unscrew and remove the two captive fixing screws.
  - c. Keeping the dummy module flap cover open, lift up the module from its seat in the injector/detector housing. Place the dummy module on a clean surface.
  - d. **DO NOT REMOVE** the gas block plug from the gas connections. See [Figure 512](#). The gas supply is done through the GDI electrical interface.

**Figure 512.** Detector Gas Block Plug (1)

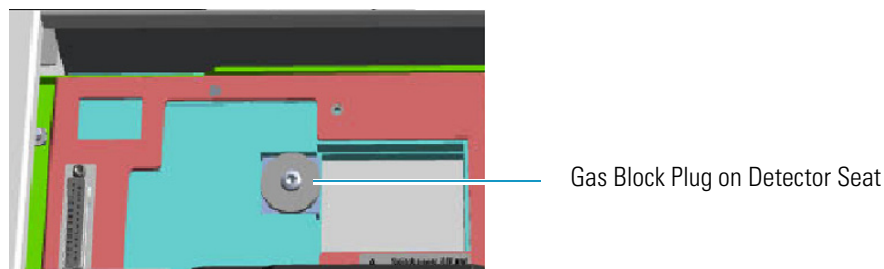

**CAUTION** In the case you are removing a real detector module instead of a dummy module, you must place and fix the gas block plug on the gas connections using a T20 Torxhead screwdriver. See [Figure 513](#).

**Figure 513.** Detector Gas Block Plug (2)

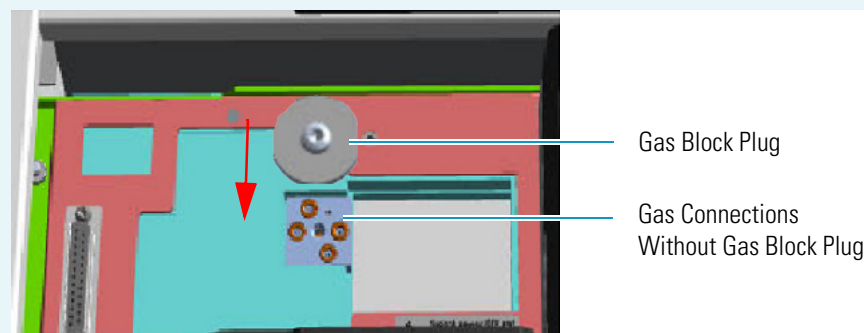

2. Plug the GDI mechanical module into the main frame.
  - a. Open the module flap cover.

## 9 Adding Modules

### Adding a Generic Detector Interface

- b. Keeping the module flap cover open, place the module in its seat. Note that the 25-pin female connector on the detector seat of the detector housing is not used.
- c. Use a T20 Torxhead screwdriver to tighten the three captive fixing screws without overtightening.



**ATTENTION** To maintain the correct alignment the screws must be tightened in turn. Tighten each screw only a small amount before moving to the next screw. Repeat until all are secure.

- d. Close the module flap cover.

## Connecting the Detector Gas Tubings to the Manifolds


This sections provides instruction for connecting the detector gas tubings to the manifolds of the GDI mechanical module and GDI electrical interface.

### ❖ To connect the gas tubing block to the manifolds

Figure 514 shows manifolds of the GDI mechanical module and GDI electrical interface.

**Note** Figure 514, Figure 515, Figure 516, Figure 517, and Figure 522 do not show the third-party detector for graphic convenience. Suppose that it is installed in the detector module GDI.


**Figure 514.** Detector Gas Tubings Manifolds (1)

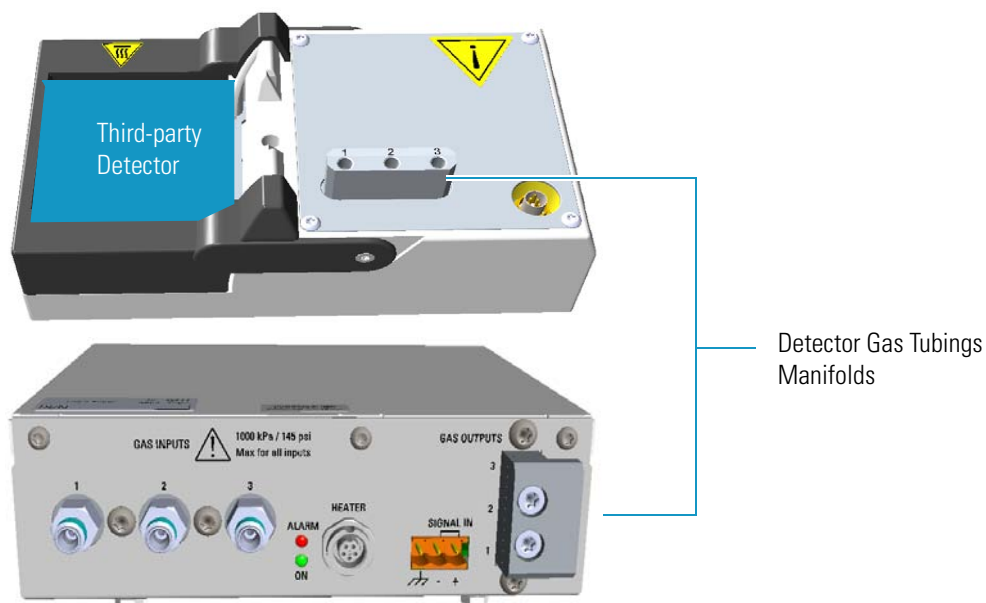

1. Take one, two, or three segments of the 1/16-in.stainless steel tubings (provided) according to the detector gases required and long enough to properly connect both the detector gas tubings manifolds.

2. Connect the detector gas tubings to the manifold of the GDI mechanical module.  
See [Figure 515](#).

**Figure 515.** Detector Gas Tubings Connections

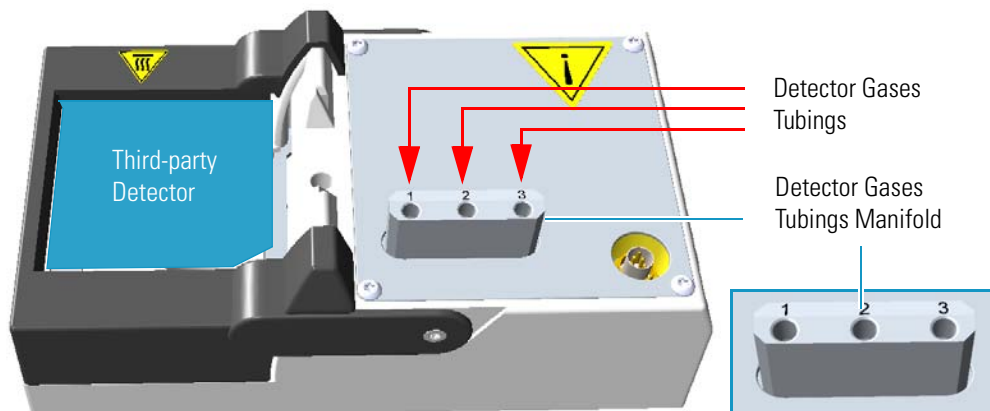

- a. Connect the detector gas tubings to the numbered inlet ports using the appropriate Swagelok® 1/16-in. nut and ferrules. Use a 1/4-in. wrench to tighten the fittings.

**ATTENTION** The inlet ports of the manifold on the GDI mechanical module are numbered **1**, **2**, and **3** respectively. Pay attention to the correct order when you connect each detector gas tubing to the corresponding outlet ports marked **Gas 1**, **Gas 2**, and **Gas 3** on the GDI electrical interface.




Connect the detector gas tubings for the third-party detector as follows:

- Outlet port **Gas 1** to inlet port 1
- Outlet port **Gas 2** to inlet port 2
- Outlet port **Gas 3** = Air (Wall) to inlet port 3


- b. Repeat [step a](#) until all the required detector gas tubings are connected to the GDI mechanical module.
- c. Bend and run the detector gas tubings along the top cover.

**Note** The bending of the detector gas tubings shown in [Figure 516](#) and [Figure 517](#) is indicative.

**Figure 516.** Bending of the Detector Gas Tubings on the TRACE 1300/1310 GC

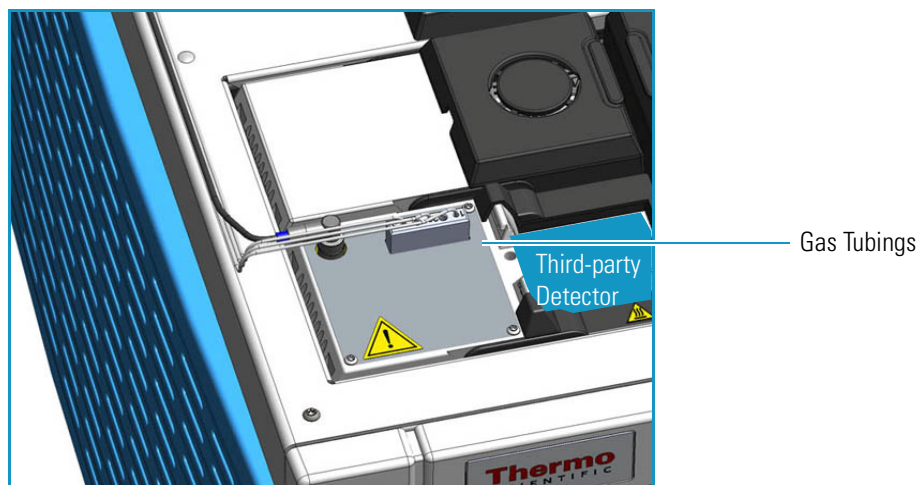


**Figure 517.** Bending of the Detector Gas Tubings on the TRACE 1310 Auxiliary Oven

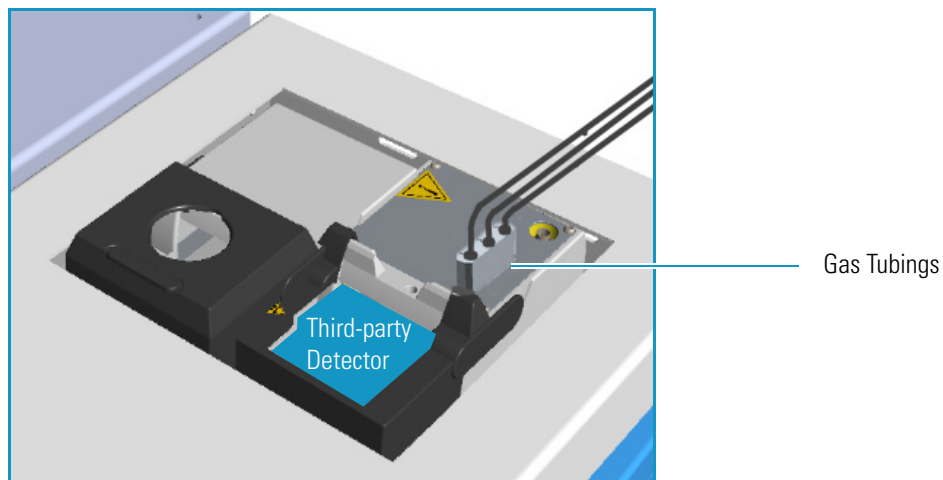

3. Connect the relevant detector gas tubings to the manifold of the GDI electrical interface. See [Figure 518](#).

**Figure 518.** Detector Gas Tubings Connection to the GDI Electrical Interface - External View (1)

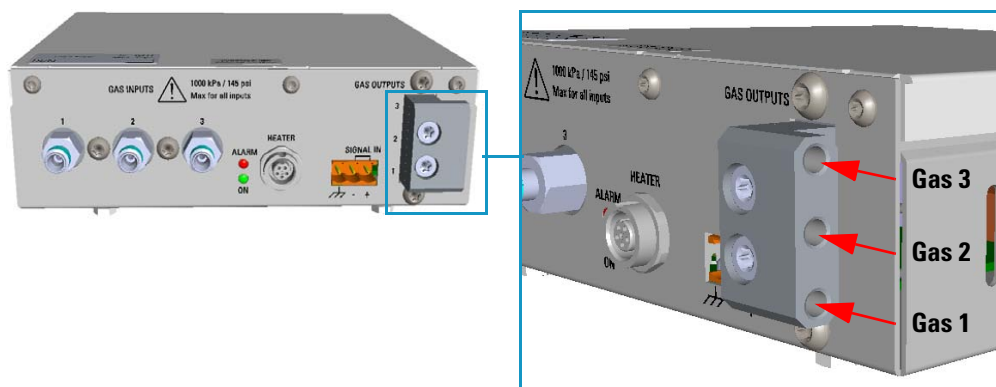


- a. Guide the detector gas tubings up to reach the GDI electrical interface.

- b. Bend and run the detector gas tubing along the back panel of the GC (or TRACE 1310 Auxiliary Oven) until its end reaches the corresponding numbered gas outlet port on the GDI electrical interface.
- c. Connect the detector gas tubing to the corresponding numbered outlet port using the appropriate Swagelok® 1/16-in. nut and ferrules. Use a 1/4-in. wrench to tighten the fittings.
- d. Repeat [step b](#) and [step c](#) until all the required detector gas tubings are connected to the GDI electrical interface.

The result of the operation is shown in [Figure 519](#).

**Figure 519.** Detector Gas Tubings Connection to the GDI Electrical Interface - External View (2)

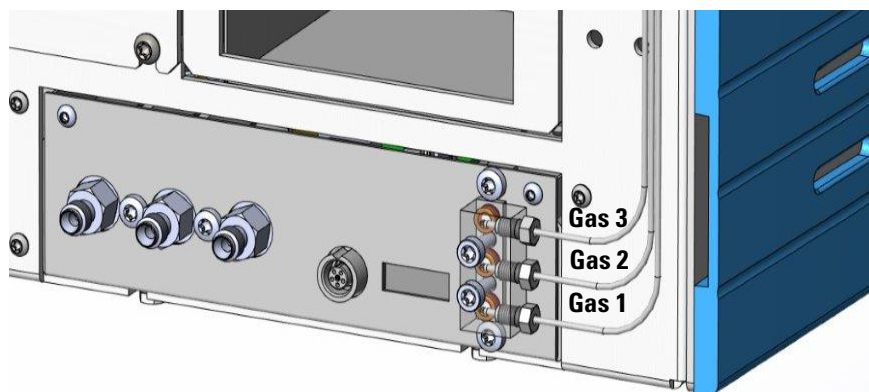


**Note** The bending of the gas tubing is indicative.

4. Connect the supply gas to the GDI electrical interface.

**Note** Use the 1/8-in. Swagelok fittings provided on the gas inlet ports to connect the gas lines.

- a. Connect the gas line to the corresponding inlet port using the appropriate nut and ferrules. Use a 7/16-in. wrench for tightening the fittings.
- b. Repeat [step a](#) until all the gas lines are connected to the corresponding inlet port on the GDI electrical interface. See [Figure 520](#).

**Figure 520.** Gas Lines Connection to the GDI Electrical Interface

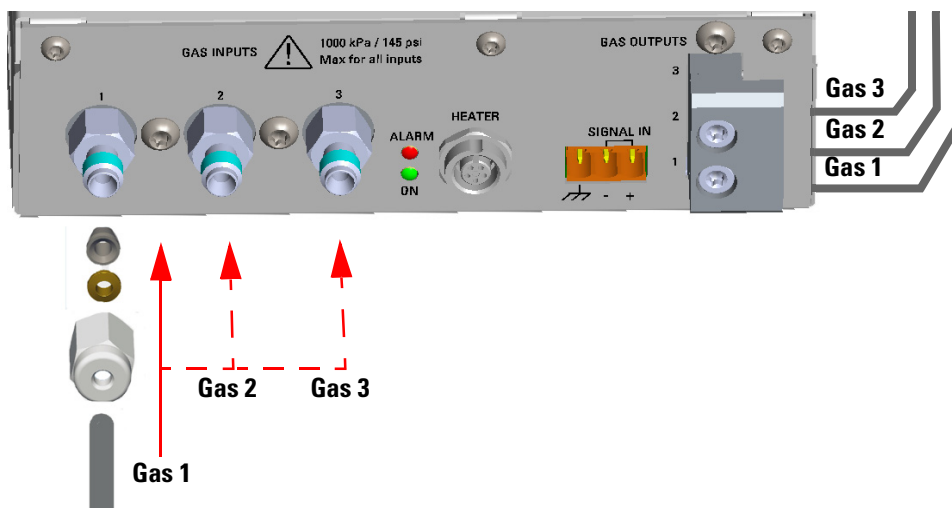

**IMPORTANT** The maximum nominal inlet pressure for all the inputs is 1000 kPa (145 psig). The working inlet pressure range is from 400 kPa (58 psig) to 1000 kPa (145 psig).


## Connecting Heater and Signal Cables

If required by the third-party detector, connect the heater cable, the signal cable, or both.


### ❖ To connect the heater and signal cables

1. Using the cable provided connect the 5-pin connector marked **Heater** on the GDI electrical interface to the detector connector on the GDI mechanical module.  
See [Figure 521](#) and [Figure 522](#).

**Figure 521.** Heater and Signal Connectors

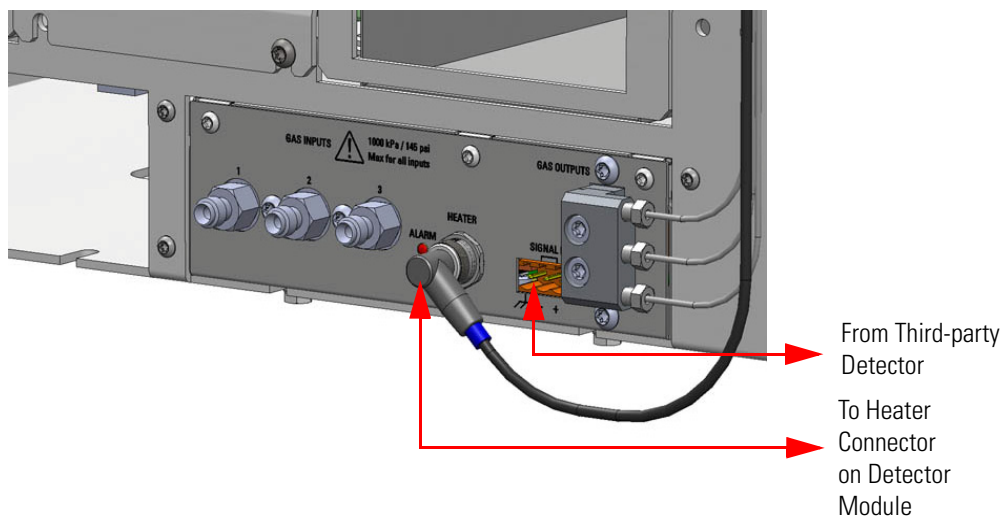


**Figure 522.** Heater Cable Connection

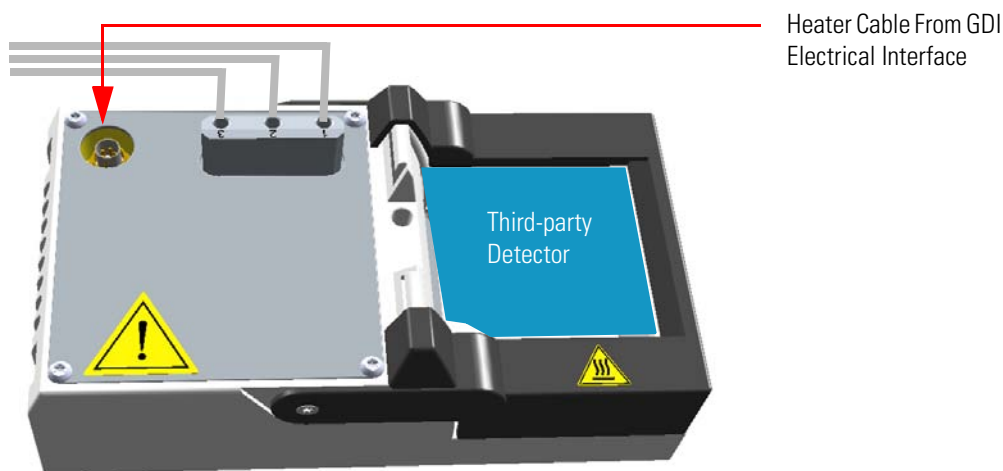

2. Connect the signals cable from the third-party detector to the connector marked Signal IN on the GDI electrical interface. See [Figure 523](#).

**Figure 523.** Signal Connection

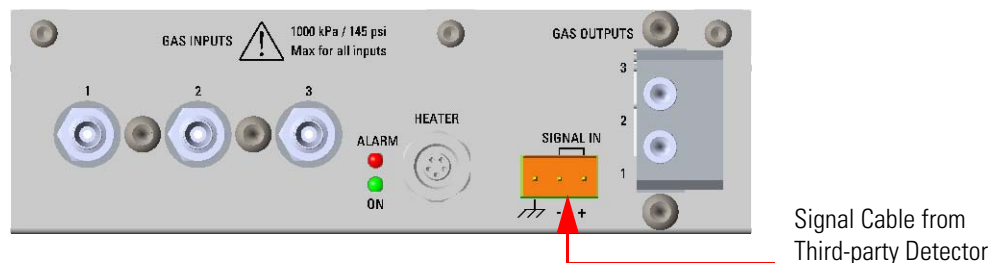

## Restarting the GC

### ❖ To restart the GC

1. Open the gas supplies.
2. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
3. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker), located at the back of the instrument, to the position I.

## Performing the Third-party Detector Start-up and Optimization

Refer to the third-party detector manual.

## Configuring and Setting GDI Detector

Configure and enable the GDI mechanical module through the user interface of your GC, or through the CDS in use. Refer to the *TRACE 1300 and TRACE 1310 User Guide*.

See the following sequences:

- “Configure and set GDI Detector through the touch screen” on page 490
- “Configure and set GDI Detector through the chromatography data system (CDS)” on page 490

### ❖ Configure and set GDI Detector through the touch screen

1. In the main menu select the Configuration icon, the Configuration menu appears. Select the Front/Back/Aux L/Aux R **GDI** detector icon to open the relevant sub-menu.
2. In the main menu select the **Instrument control** icon. The Instrument Control menu appears. In the Instrument Control menu, select the Front/Back/Aux L/Aux R GDI detector icon to open the relevant sub-menu.
3. Set the **parameters** values as required, then return to main menu.

**Note** For details refer to **Chapter 2** in the *TRACE 1300/REACE 1310 User Guide*.

### ❖ Configure and set GDI Detector through the chromatography data system (CDS)

1. In the **Configuration** window select the **Detectors** tab.
2. Select the **Detector Type**: choose **GDI**.
3. Click **GDI Config...**; the GDI Configuration page is visualized.
4. Select the **Gas Type** used for the detector gases **Gas 1**, **Gas 2**, and **Gas 3**. Choose one: Air, Hydrogen, Nitrogen, Helium, Argon, or Argon/Methane. Nitrogen is the default gas.
5. Set the **Full-scale flow** of the restriction installed for each detector gas. Set a value in the range 1-1000 mL/min. Default value is 50 mL/min.
6. Select the **Max. detector temperature** in the range from 0 °C to 450 °C. The default temperature is 400 °C.
7. Select the **ADC full-scale voltage**. Choose one: 1 V, 5 V, or 10 V. The default value is 1V.
8. Open the **GDI Method Page** and set the required parameters.

**Note** For details refer to **Chapter 4** in the *TRACE 1300/REACE 1310 User Guide*.

**Note** If the heater is installed but the actual Temperature read back is 0 °C, check the integrity of the heater and the temperature probe.



## Adding an Analog Output Interface

This section provides instructions for installing the Analog Output Interface (AOI) on your TRACE 1300/1310 GC.



**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.


### ❖ To install an analog output interface

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


3. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
4. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
1. Remove the cover of the external modules housing where installing the Analog Output Interface module. See [Figure 524](#).

**Figure 524.** Housing Cover Removal

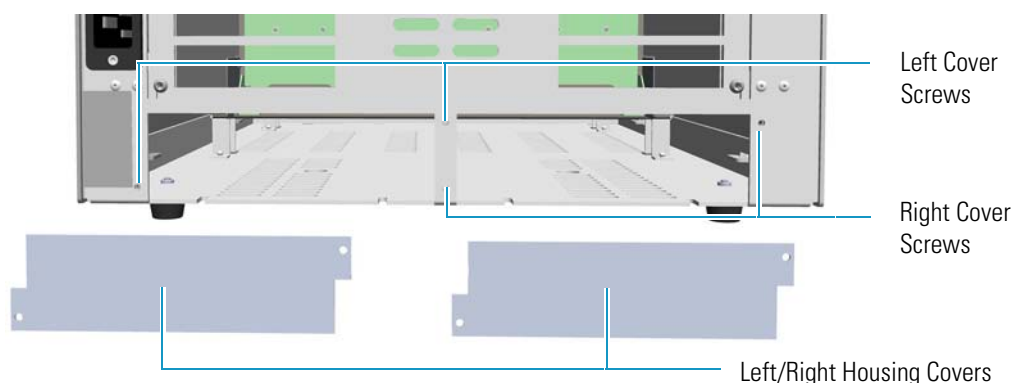

- a. Using a T20 Torxhead screwdriver, unscrew and remove the left or right housing cover screws, then remove the cover from the housing.
2. Install the Analog Output Interface into the housing.
  - a. Loosen the two hexagonal screws under the module. See [Figure 525](#).


**Figure 525.** AOI Installation (1)

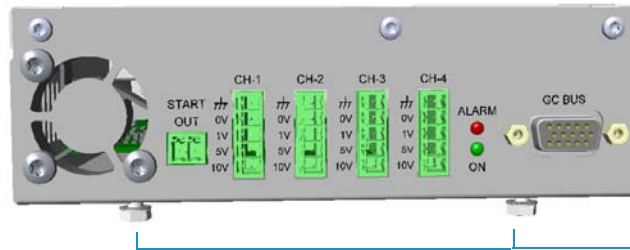

Hexagonal Screws

- b. Carefully place the Analog Output Interface into the left or right housing until the hexagonal screws couple with the slots on the floor of the GC. See [Figure 526](#).

**Figure 526.** AOI Installation (2)

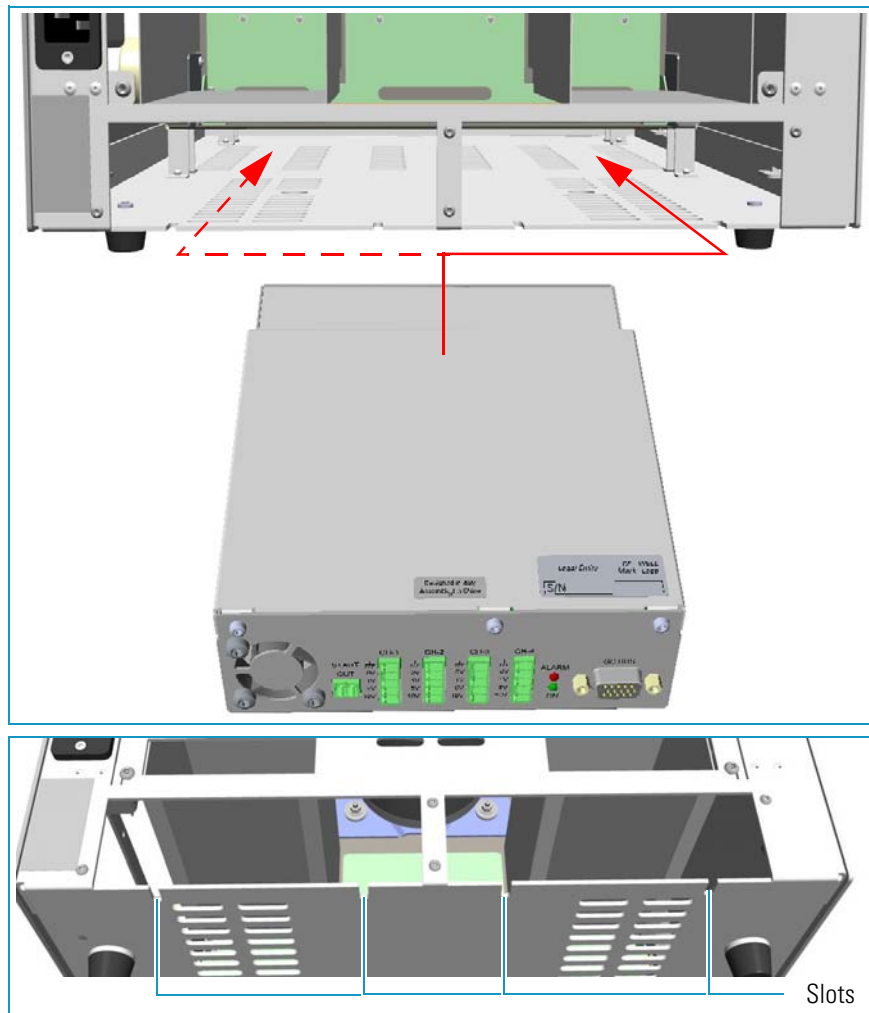


Slots

- c. Finger-tighten the hexagonal screws slightly, or use a 10-mm wrench.

**Note** Always keep the hexagonal screws in their place. This allows you an easier removal of the auxiliary module when necessary.


The result of the installation is shown in [Figure 508](#).

**Figure 527.** Analog Output Interface Module Installed into the GC

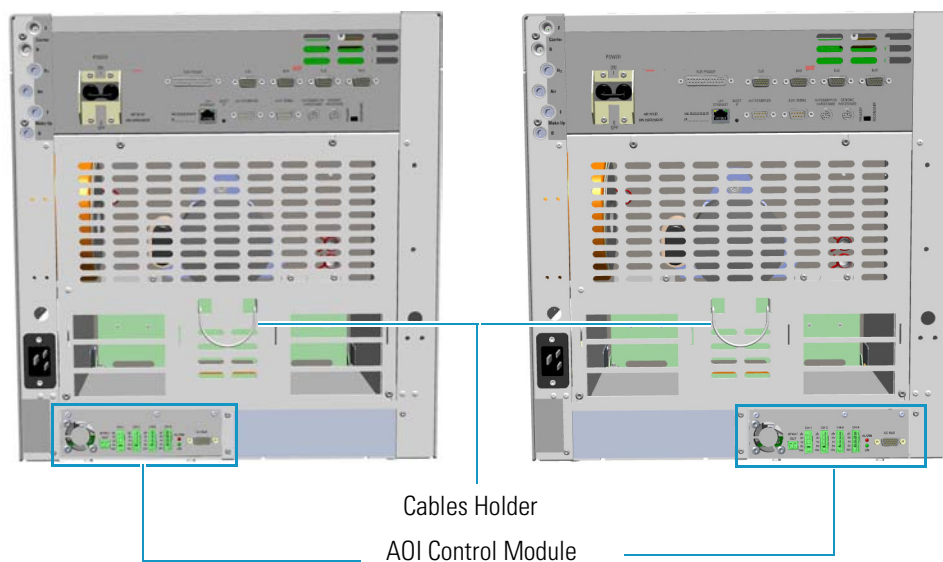

3. Connect the analog output interface module. See [Figure 528](#).

**Figure 528.** AOI Connections

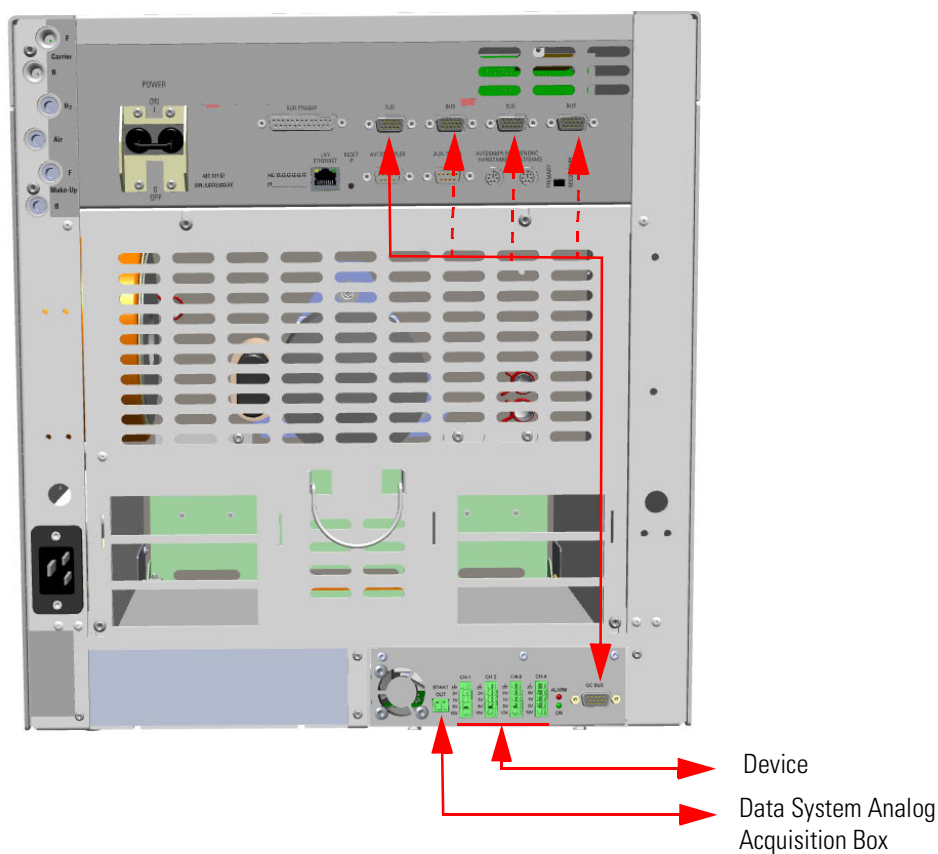

**Tip** To avoid the contact with the hot air from the vents, it is suggested to gather the electrical cables into cables holder.

- a. Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.
- b. Connect up to four channels **CH-1**, **CH-2**, **CH-3**, **CH-4** to your device selecting the full scale of each analogue signal from **1 V** to **10 V** according to your needs.



**IMPORTANT** Only connect the Ground wire to the GND contact of the AOI module or of the your device. Do not connect the Ground cable to both the devices.

- c. Connect the **Start OUT** contact closure to **Start IN** TTL line of the Data system analog acquisition box.
4. Open the gas supplies.
  5. If other external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
  6. Power on the GC.
    - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
    - b. Flip up the power switch (breaker), located at the back of the instrument, to the position I.
  7. Configure and enable the system installed through the touch screen of your GC. Refer to the *TRACE 1300 and TRACE 1310 User Guide*.

## Adding Systems

This chapter describes how to install any added system that is available for the TRACE 1300/TRACE 1310. See the *TRACE 1300 and TRACE 1310 Spare Parts Guide* for information about ordering the equipment in this chapter.


### Contents

- [Adding the Oven Cryo System](#)
- [Adding the PTV and PTVBKF Cryo System](#)
- [Adding an Auxiliary Gas System](#)
- [Adding the Hydrogen Sensor](#)


## Adding the Oven Cryo System

This section provides instructions for installing and configuring the Oven Cryo system on your TRACE 1300/TRACE 1310 using the dedicated kit. See [Figure 529](#).

**Figure 529.** Oven Cryo System for Carbon Dioxide and Liquid Nitrogen

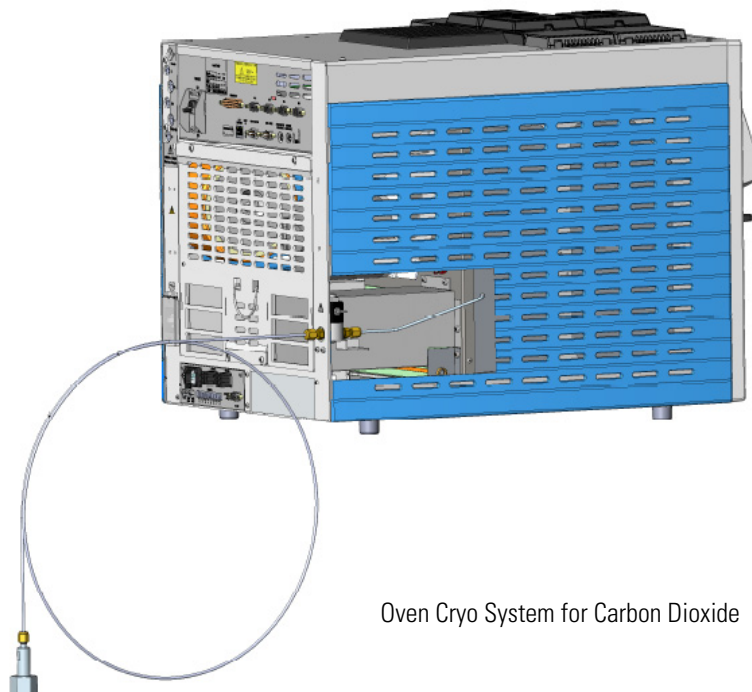

Oven Cryo System for Carbon Dioxide

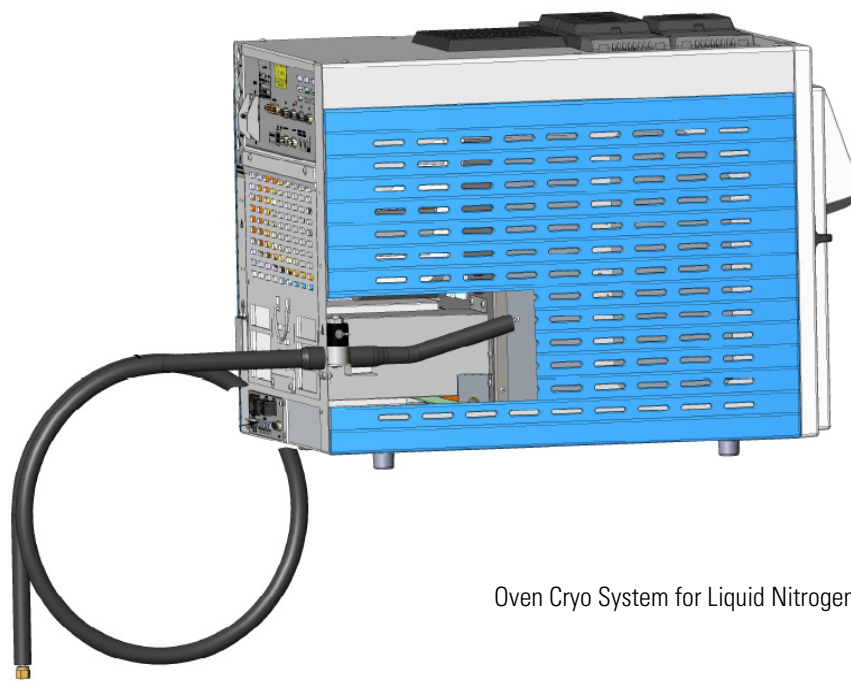

Oven Cryo System for Liquid Nitrogen

**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.



The cryo system requires the use of liquid nitrogen or carbon dioxide as coolant. Before using liquid nitrogen or carbon dioxide, read the indications of hazards, and the instructions reported in the Safety sheet supplied by the manufacturer with reference to the CAS number (Chemical Abstract Service). See also “[Liquid Nitrogen Safety Precautions](#)” on [page xxix](#) and “[Carbon Dioxide Safety Precautions](#)” on [page xxx](#).


## Oven Cryo System Overview

Two Oven Cryo Upgrade Kits are available:


- Oven Cryo with Liquid Nitrogen (LN<sub>2</sub>) as coolant.
- Oven Cryo with Carbon Dioxide (CO<sub>2</sub>) as coolant.

Each upgrade kit contains all the material required to install the Oven Cryo system on your GC. See [Figure 530](#) and [Figure 531](#).

**Figure 530.** Oven Cryo Kit for Carbon Dioxide

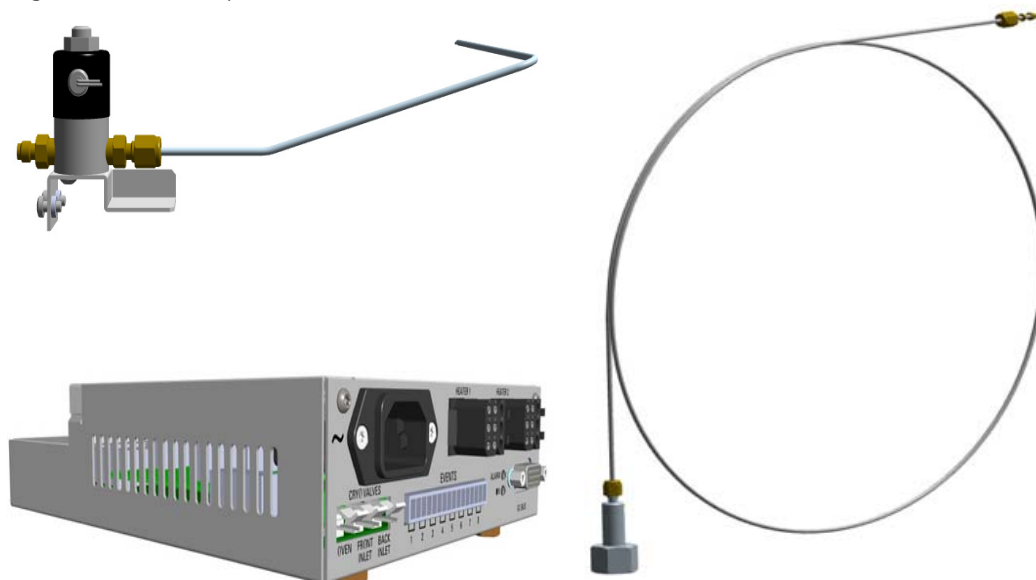


**Figure 531.** Oven Cryo Kit for Liquid Nitrogen

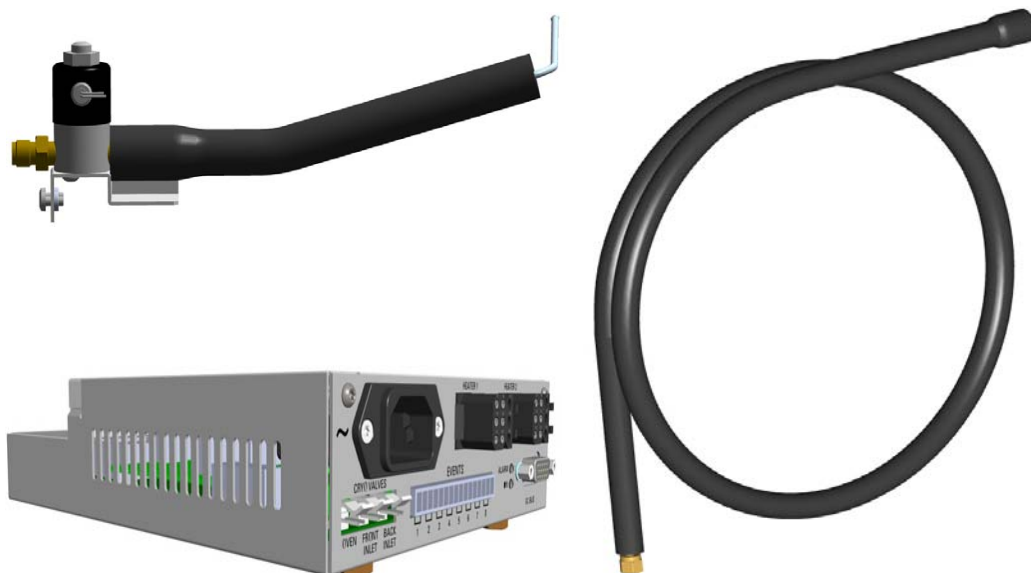

Each upgrade kit contains all the material required to install the Oven Cryo system on your GC:

- Dedicated solenoid valve mounted on a support bracket.
- Tube for the coolant into the oven.
- Coolant tank tube with connection fittings.
- Aux Temperature/Cryo Module that should be installed into a free external module housing provided on the back of the GC. See [Figure 532](#).

**Figure 532.** External Modules Housing

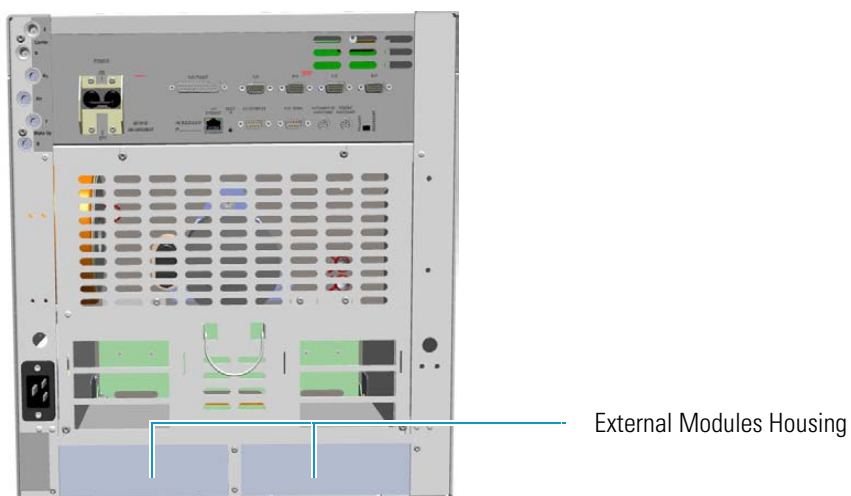



## Installing the Oven Cryo System


### ❖ To install the Oven Cryo System

1. Put the GC in standby condition.
2. Cool the oven, injectors, and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.


3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, in the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
6. Remove the Left Side panel.
  - a. Open the front door of the GC. Using a T20 Torxhead screwdriver, unscrew the left side panel screw from the interior front panel. See [Figure 533](#).

**Figure 533.** Left Side Panel Fixing Screw

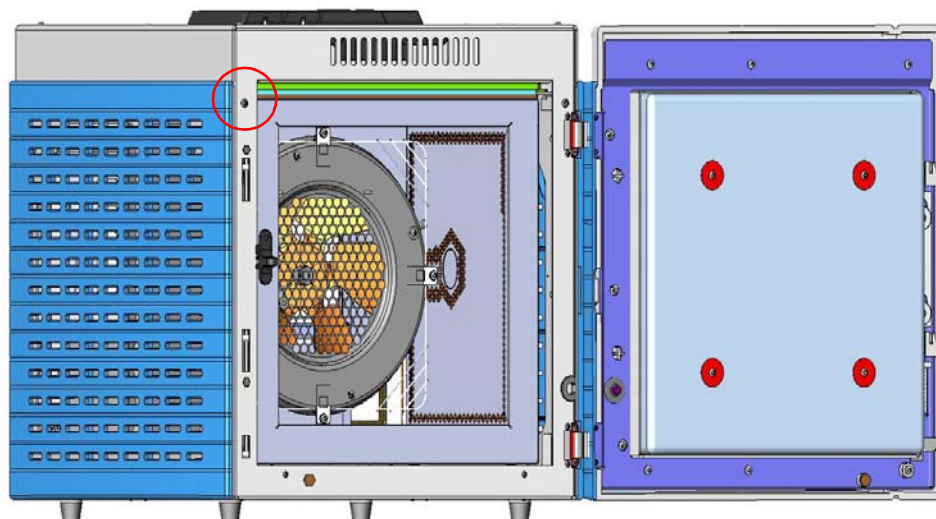

- b. Slide the panel towards the back of the instrument up to the stop.
- c. Remove the panel by pulling it outwards. Be aware that the ground wire is attached to the panel. See [Figure 534](#).

**Figure 534.** Left Panel Removal

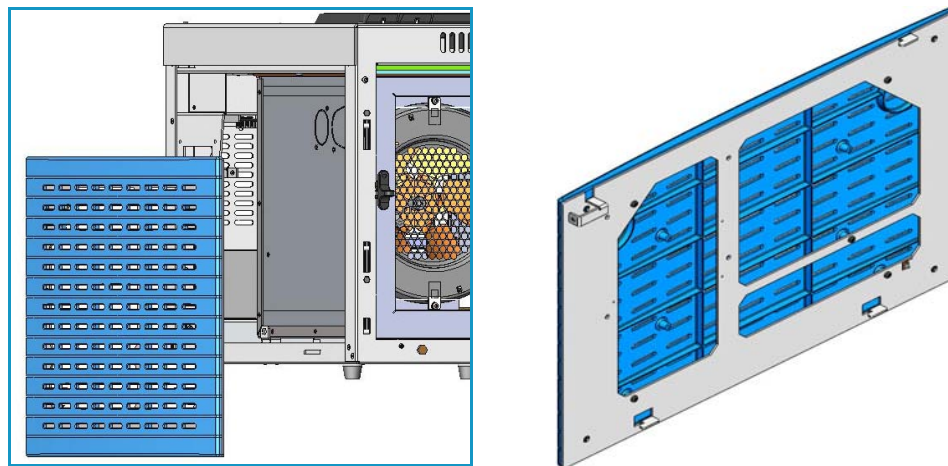


- d. Unplug the ground wire from the panel.

**Note** Pay attention to the positioning of the ground wire plug, so it can be reconnected in the same way it was removed.


7. Install the Aux Temperature/Cryo Module. See [Figure 535](#)

**Figure 535.** Housing Cover Removal

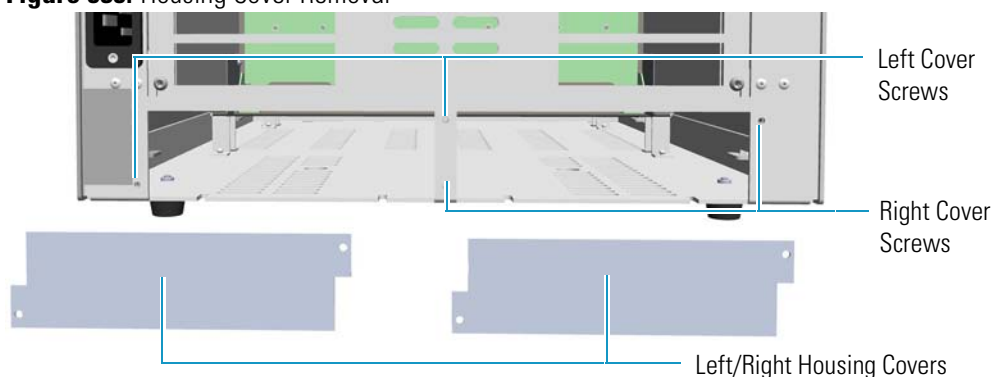


- a. Remove the cover of the external modules housing where installing the module.
  - b. Using a T20 Torxhead screwdriver, unscrew and remove the left or right housing cover screws.
  - c. Remove the covers from the housing.
8. Install the module into the housing
    - a. Loosen the two hexagonal screws under the module. See [Figure 536](#).

**Figure 536.** Module Installation (1)

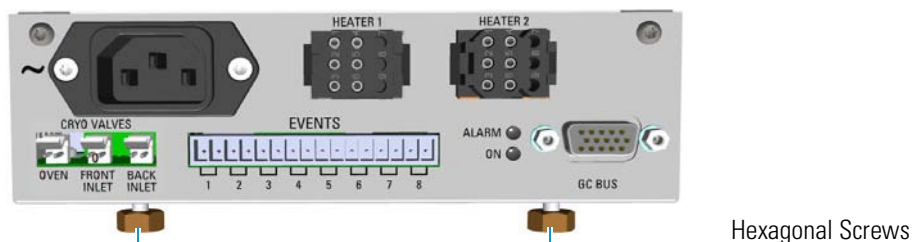

- b. Place the module into the left or right housing until the hexagonal screws couple the slots on the floor of the GC. See [Figure 537](#).

**Figure 537.** Module Installation (2)

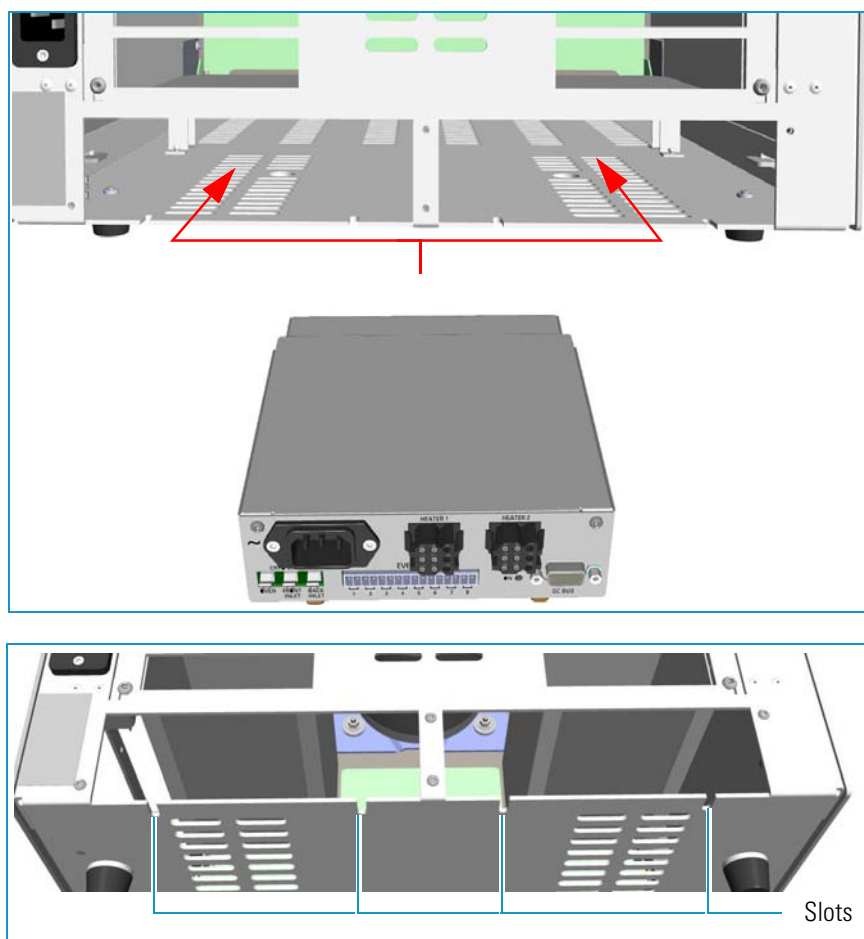


- c. Finger-tighten the hexagonal screws slightly, or use a 10-mm wrench.

**Note** Always keep the hexagonal screws in their place. This allows you an easier removal of the auxiliary module when necessary.

The result of the installation is shown in [Figure 538](#).


**Figure 538.** Aux Temperature/Cryo Module Installed into the GC

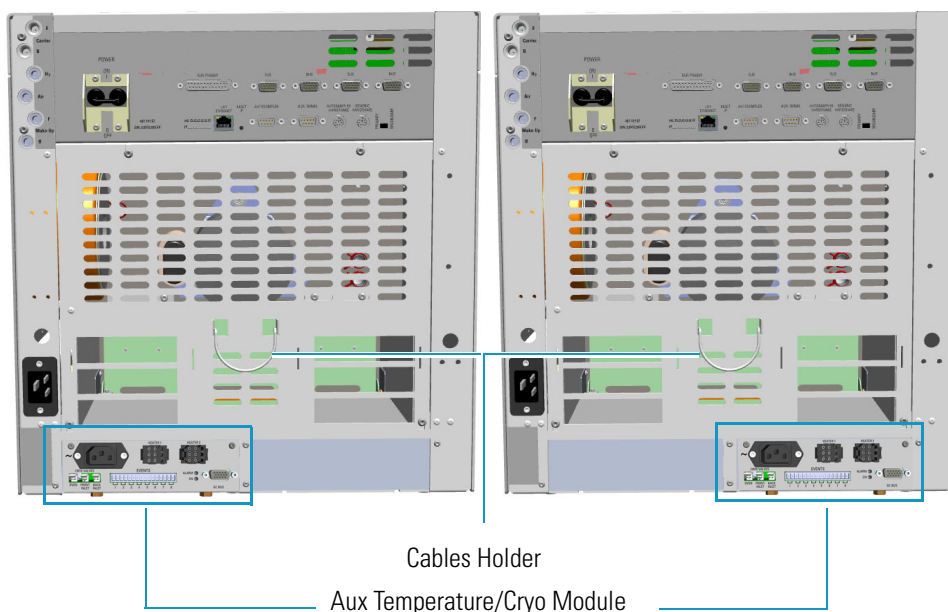

9. Install the solenoid valve.

The solenoid valve must be installed in the proper seat on the back of the GC. See [Figure 539](#).


**Figure 539.** Solenoid Valve Assembly Seat

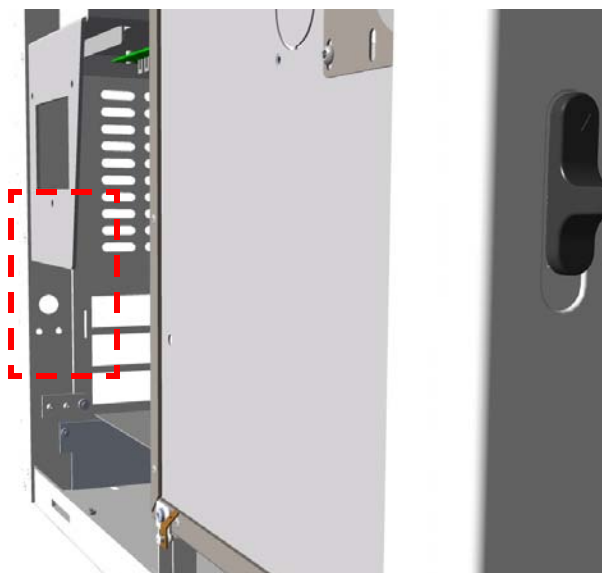

- a. With care, take the solenoid valve assembly and remove the two screws from the brackets using a T15 Torxhead screwdriver. These screws are used to fix the bracket on the GC.


If not already done, connect the coolant tube to the solenoid valve using the proper nut and ferrule. Use a 7/16-in. wrench for tightening the fittings. See [Figure 540](#).

**Figure 540.** Solenoid Valve Assembly

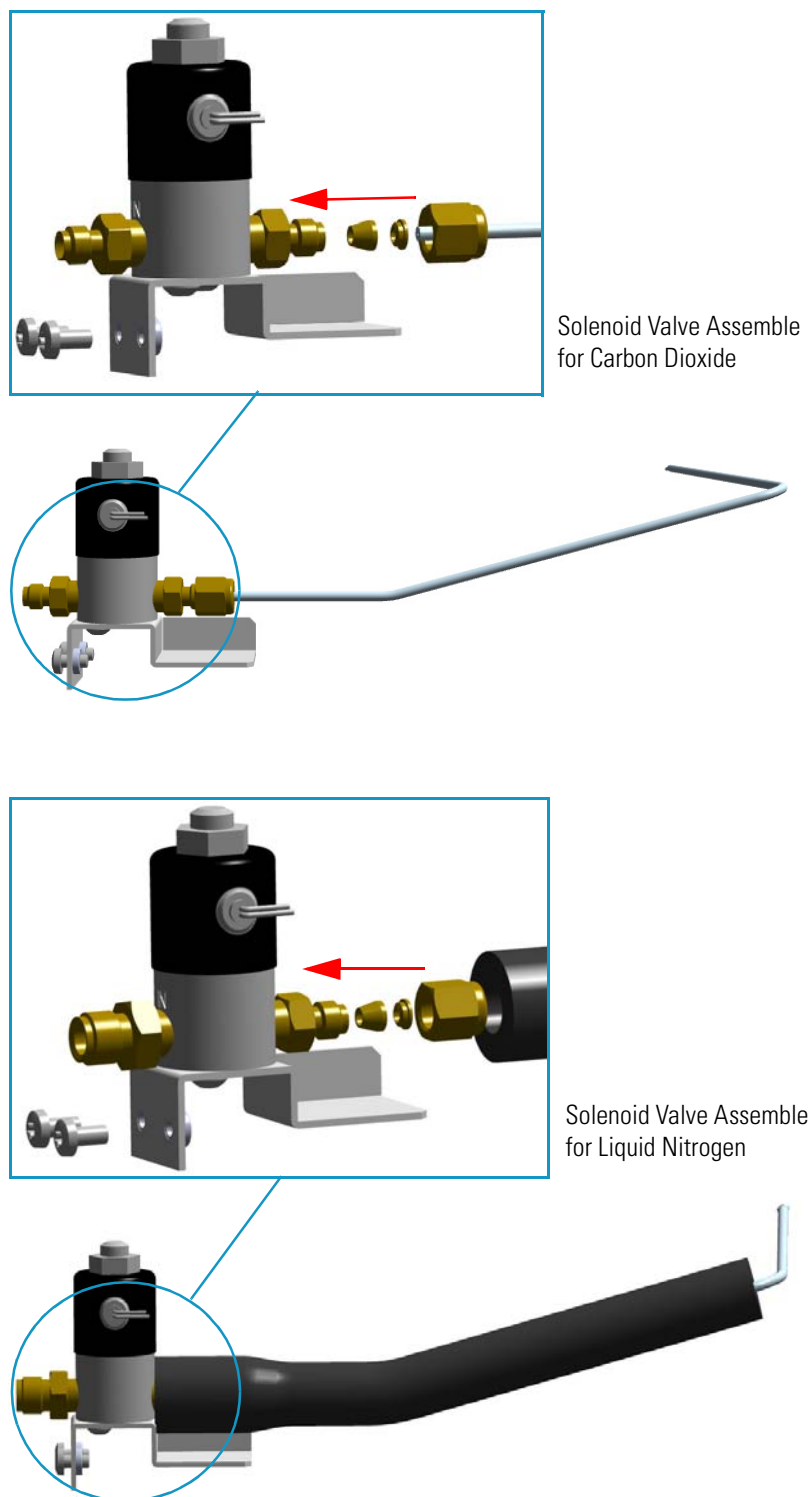


- b. Guide the solenoid valve assembly into its seat on the back of the GC, and the tube for the coolant into the oven through the holes provided. See [Figure 541](#) and [Figure 542](#).

**Figure 541.** Installation of the Solenoid Valve Assembly for Carbon Dioxide

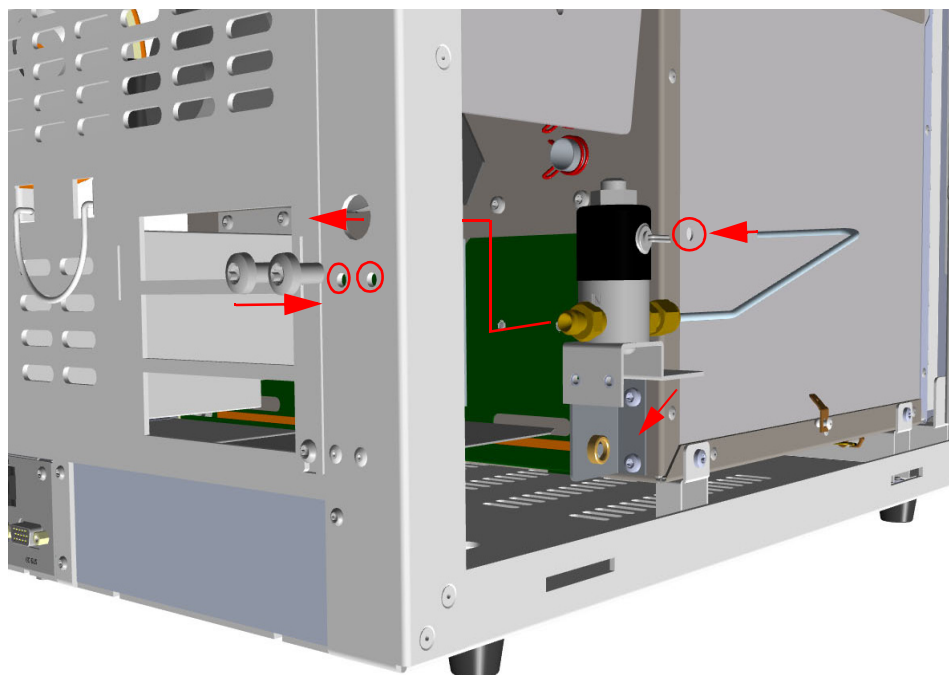

**Figure 542.** Installation of Solenoid Valve Assembly for Liquid Nitrogen

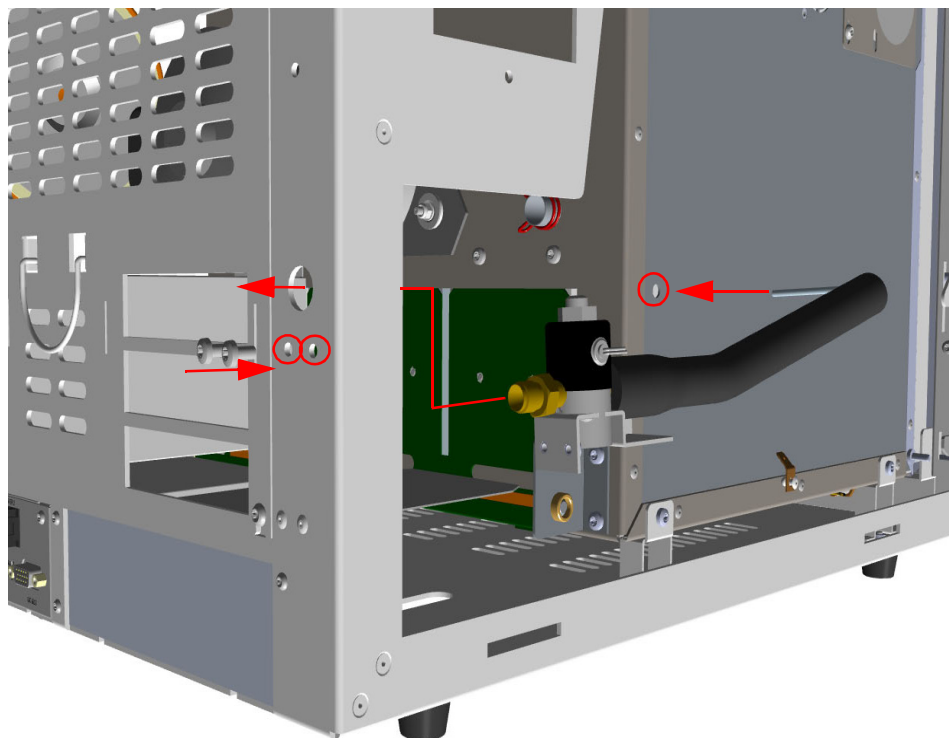

- c. Fix the solenoid valve assembly to the back of the GC using the two fixing screws previously removed. See [Figure 543](#).


**Figure 543.** Solenoid Valve Assemble Fixing



10. Connect the cryogenic tank tube to the solenoid valve assembly.
  - a. Connect the proper end of the cryo supply tube to the 1/8-in. NPT connection of the solenoid valve using the proper nut and ferrule. Use a 7/16-in. wrench for tightening the fittings. See [Figure 544](#) and [Figure 545](#).




**Figure 544.** Oven Cryo Supply Tube for Carbon Dioxide Connection

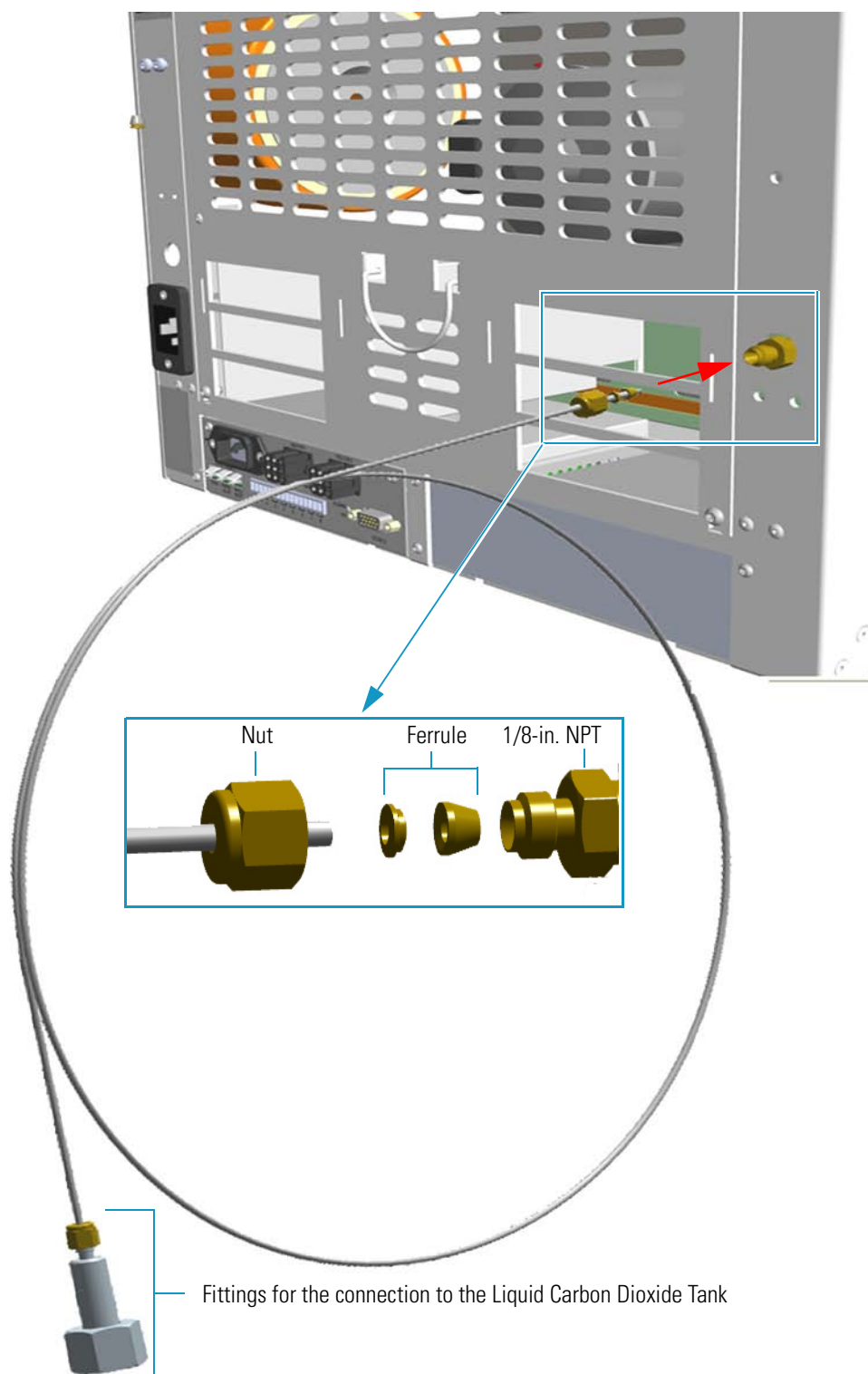




**Figure 545.** Oven Cryo Supply Tube for Liquid Nitrogen Connection

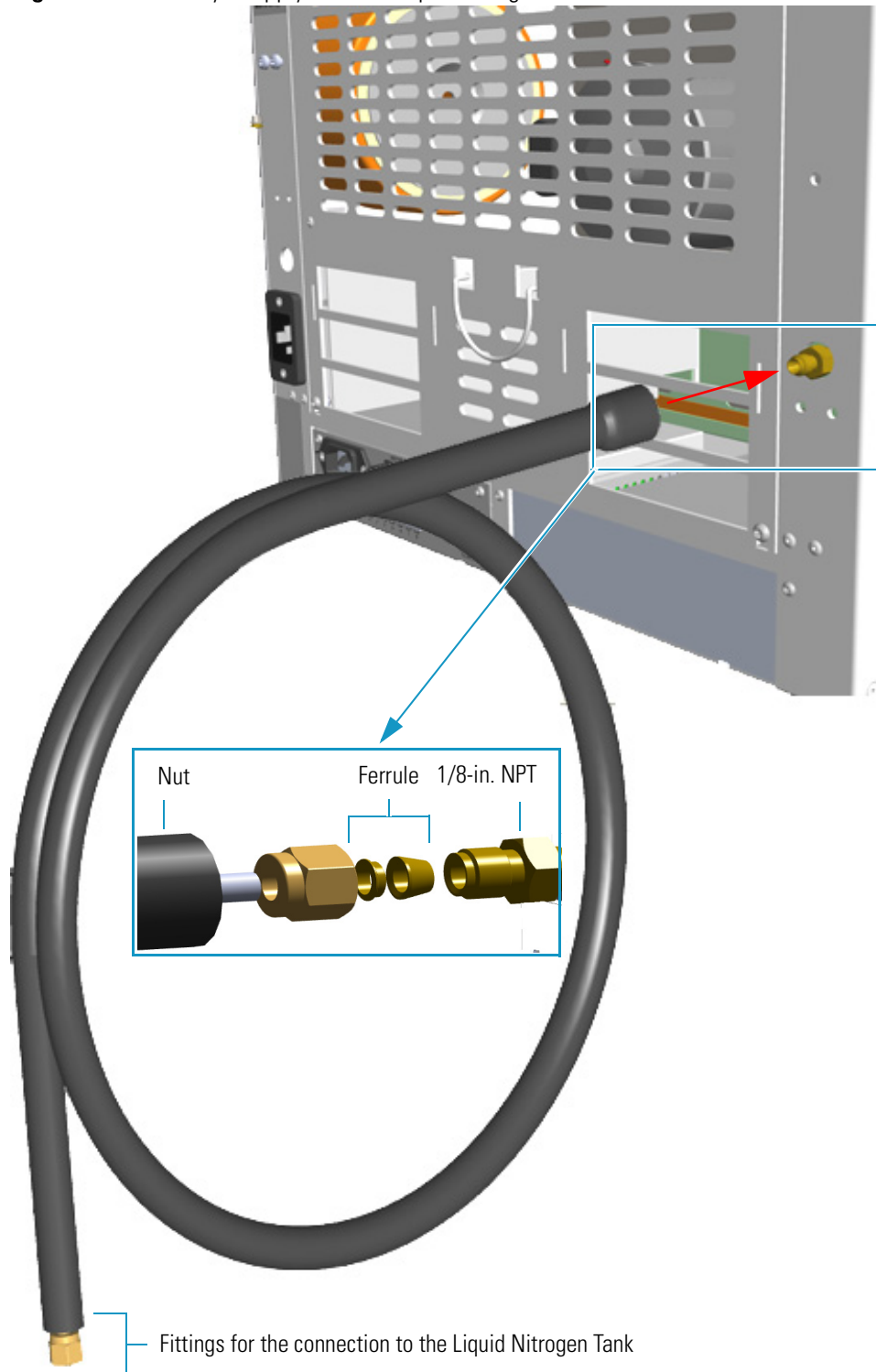

- b. Connect the other end to the coolant container using the appropriate nuts and ferrules.


11. Connect the cryo valve to the Aux Temperature/Cryo Module.
  - a. Connect the cable provided to the solenoid valve connector, and guide the cable through the slot in the centre of the back panel. See [Figure 546](#).

**Figure 546.** Solenoid Valve Cable Connection

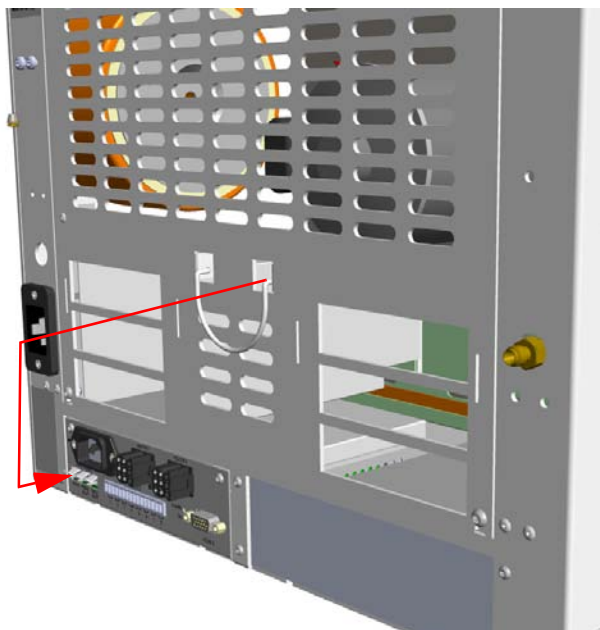

- b. Connect the cryo solenoid valve to the 2-pin connector marked **Cryo Valves - Oven** using the cable provided. See [Figure 547](#).

**Figure 547.** Cryo Valves: Oven

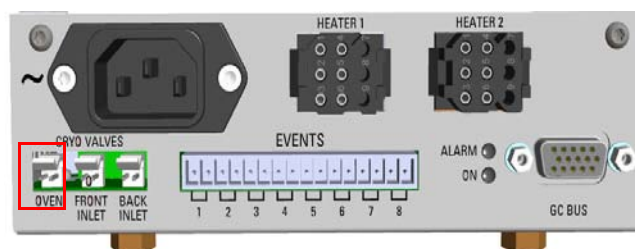

- c. Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.
  - d. Plug the power cable to the AC Input connector on the front of the module and to the wall outlet. The LED marked **On** lights after the GC is powered on.

**Note** For further details regarding the installation of the Aux Temperature/Cryo module, see the section [“Adding an Aux Temperature/Cryo Module”](#) on [page 443](#).

- e. Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.

12. Mount the left side panel proceeding in the reverse order in which the left side panel was removed.
13. Open the gas supplies.
14. If other external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
15. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker), located at the back of the instrument, to the position I.
16. Configure and enable the cryogenic system through the user interface of your GC, or through the CDS in use. Refer to the *TRACE 1300 and TRACE 1310 User Guide*.
  - a. **Configuration and enabling through the touch screen.**
    - i. In the main menu select the **Configuration** icon. The Configuration menu appears.
    - ii. In the Configuration menu, select the **Oven** icon to open the relevant submenu.
    - iii. Set the **Cryogenic** parameters.
      - **Cryogenic Type** — Enable or disable the cryogenic system when it is installed and configured with Carbon Dioxide or Liquid Nitrogen as a coolant. Select between LN<sub>2</sub>, CO<sub>2</sub>, or none.
      - **Cryo timeout** — Enter the time at which the cryo system will be disabled. Enter a value from 0–30 min.
      - **Cryo Threshold** — Specify the temperature at which the cryo system begins to supply the coolant. Enter a value from 40–200 °C.
    - iv. Return to main menu.
  - b. **Configuration and enabling through the Chromatography Data System.**
    - i. In the **Configuration** window select the **Auxiliary** tab.
    - ii. Select the **Auxiliary control module** check box to enable the setting for the auxiliary control of the module option installed on your GC.
    - iii. Select the **Oven cryogenics** check box.
    - iv. Select the **Cryo type** used by your cryogenic option. Choose one: Liquid Nitrogen or Carbon Dioxide.
    - v. Open the **Oven** page. Select the **Cryogenics enable** check box to enable the cryogenic system.
    - vi. In **Cryo threshold** text box specify the temperature at which the cryo system begins to supply the coolant. Enter a value from 40–200 °C.

## 10 Adding Systems

Adding the PTV and PTVBKF Cryo System

17. Set the normal detector, injector, and GC working conditions.


### Adding the PTV and PTVBKF Cryo System

This section provides instructions for installing and configuring the PTV and PTVBKF Cryo system on your TRACE 1300/TRACE 1310 using the dedicated kit. See [Figure 548](#) and [Figure 549](#).

**Figure 548.** PTV and PTVBKF Single and Double Cryo System for Carbon Dioxide



**Figure 549.** PTV and PTVBKF Single and Double Cryo System for Liquid Nitrogen

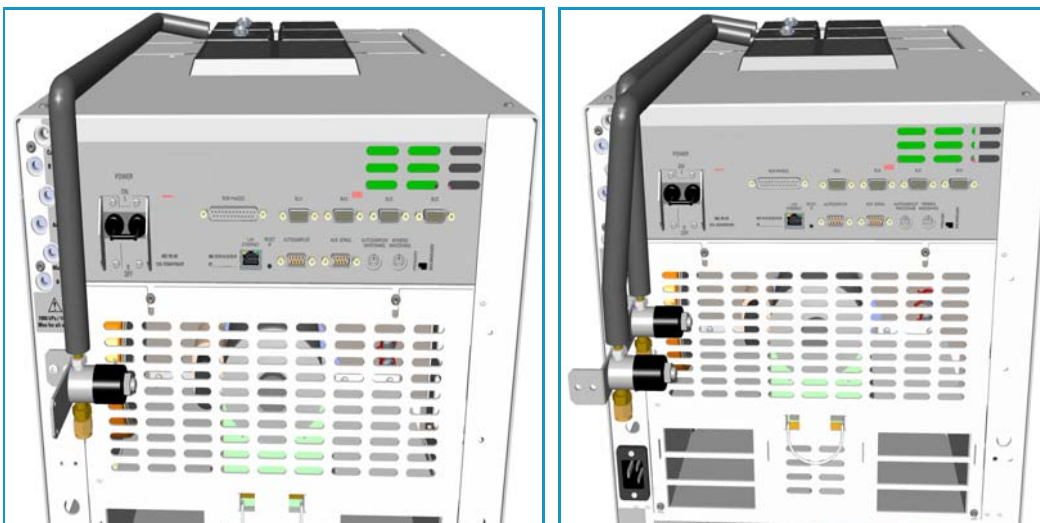



**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.

The cryo system requires the use of liquid nitrogen or carbon dioxide as coolant. Before using liquid nitrogen or carbon dioxide, read the hazard indications and the instructions reported in the Safety sheet supplied by the manufacturer with reference to the CAS number (Chemical Abstract Service). See also “[Liquid Nitrogen Safety Precautions](#)” on page xxix and “[Carbon Dioxide Safety Precautions](#)” on page xxx.

## PTV/PTVBKF Cryo System Overview

Two PTV/PTVBKF Cryo Upgrade Kits are available:

- PTV/PTVBKF Cryo with Liquid Nitrogen (LN<sub>2</sub>) as coolant.
- PTV/PTVBKF Cryo with Carbon Dioxide (CO<sub>2</sub>) as coolant.

Each upgrade kit contains all the material required to install the PTV/PTVBKF Cryo system on your GC. See [Figure 550](#) and [Figure 551](#).

**Figure 550.** PTV/PTVBKF Cryo Kit for Carbon Dioxide



## 10 Adding Systems

Adding the PTV and PTVBKF Cryo System


**Figure 551.** PTV/PTVBKF Cryo Kit for Liquid Nitrogen



Each upgrade kit contains all the material required to install the PTV/PTVBKF Cryo system on your GC:

- Dedicated solenoid valve mounted on a support bracket.
- Tube for the coolant into the PTV/PTVBKF injector module.
- Coolant tank tube with connection fittings.
- Aux Temperature/Cryo Module that should be installed into a free external module housing provided on the back of the GC. See [Figure 552](#).

**Figure 552.** External Modules Housing

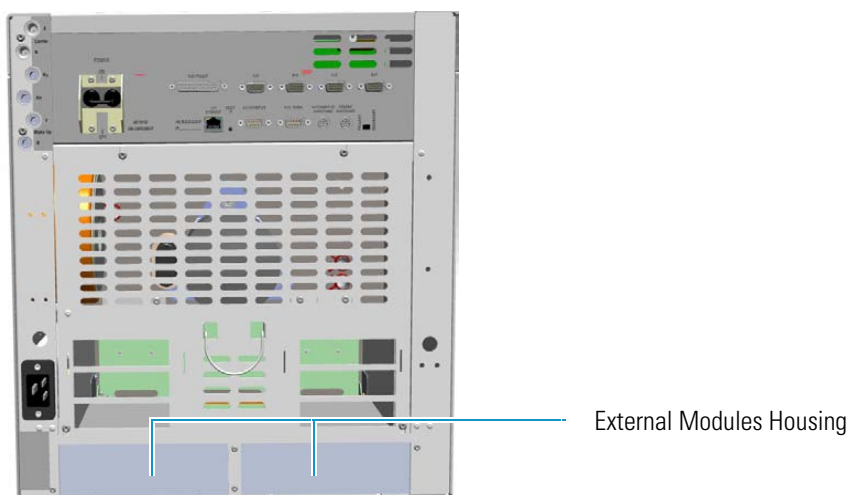

## Installing the PTV/PTVBKF Cryo System


### ❖ To install the PTV/PTVBKF Cryo System

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
6. Put the autosampler away if present.
7. Install the Aux Temperature/Cryo Module.
  - a. Remove the cover of the external modules housing where installing the module. See [Figure 553](#).


**Figure 553.** Housing Cover Removal

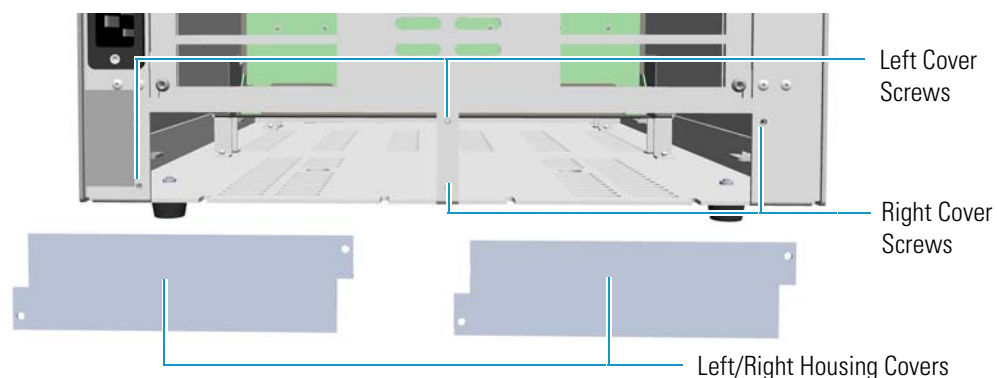

- b. Using a T20 Torxhead screwdriver, unscrew and remove the left or right housing cover screws.
  - c. Move and drive the cover out from the housing.
8. Install the module into the housing
  - a. Loosen the two hexagonal screws under the module. See [Figure 554](#).




**Figure 554.** Module Installation (1)

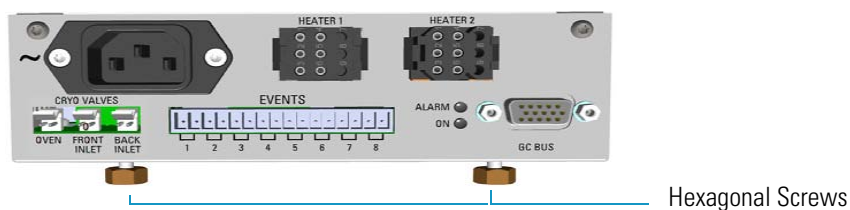

- b. Place the module into the left or right housing until the hexagonal screws couple with the slots on the floor of the GC. See [Figure 555](#).

**Figure 555.** Module Installation (2)

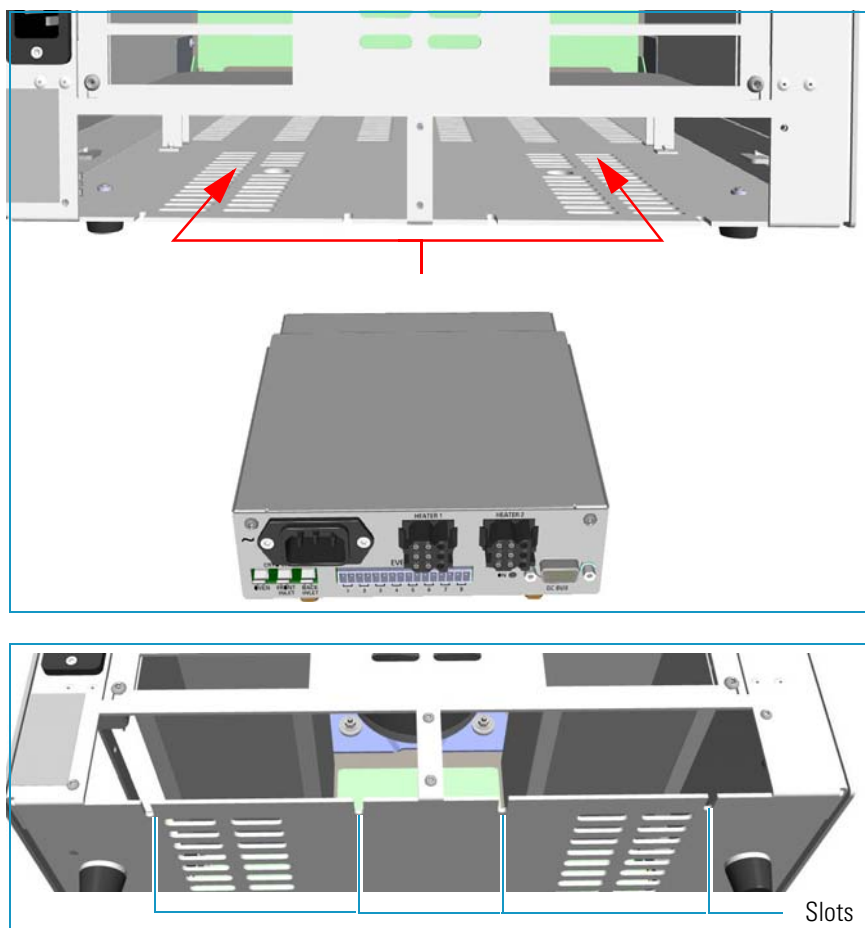


- c. Finger-tighten the hexagonal screws slightly or use a 10-mm wrench.

**Note** Always keep the hexagonal screws in their place. This allows you an easier removal of the auxiliary module when necessary.

The result of the installation is shown in [Figure 556](#).



**Figure 556.** Aux Temperature/Cryo Module Installed into the GC

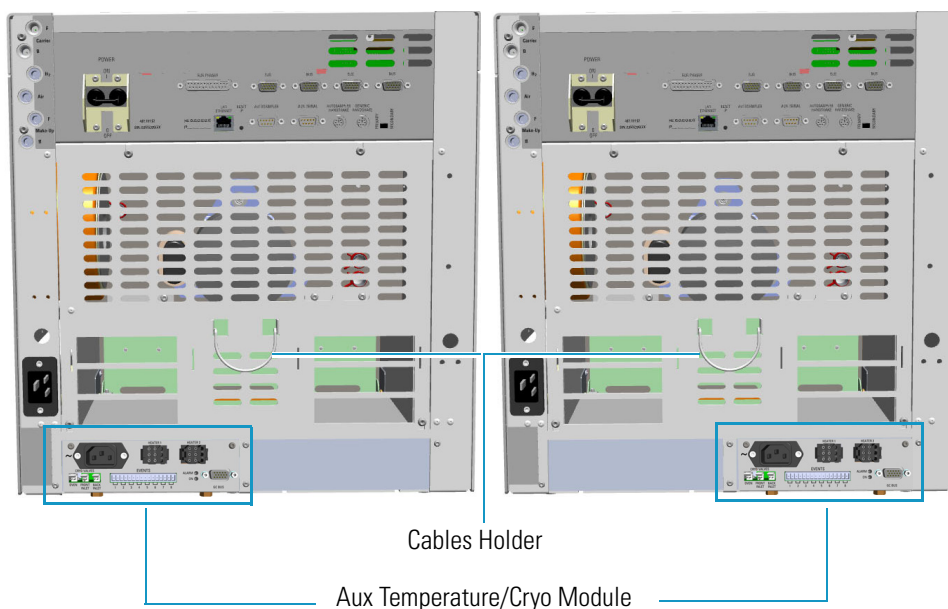


9. Install the proper solenoid valve assembly.

The solenoid valve must be installed in the proper seat on the back of the GC.

- In case of a single cryo kit, mount the solenoid assembly as shown in [Figure 557](#).
- In case of double cryo kit, mount the two solenoid valves assembly as shown in [Figure 558](#).


**Figure 557.** Single Solenoid Valve Assemble Installation

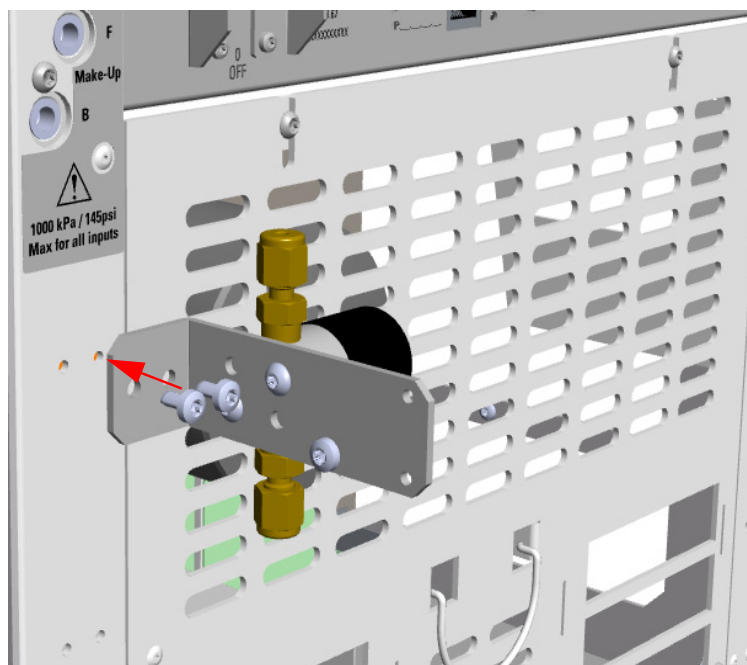

## 10 Adding Systems


Adding the PTV and PTVBKF Cryo System

**Figure 558.** Double Solenoid Valve Assemble Installation

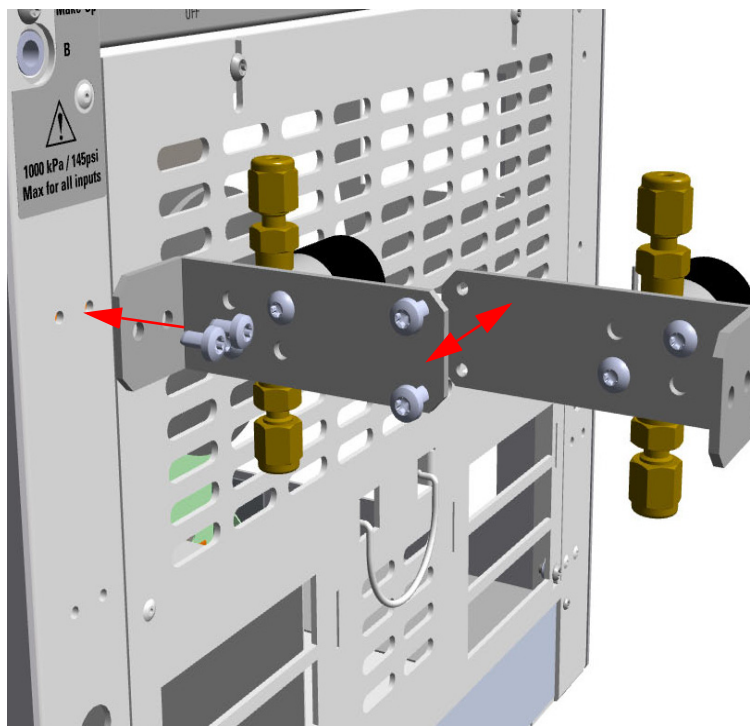


10. Insert the coolant tube into the PTV/PTVBKF injector module.
  - a. Open the module flap cover.
  - b. On the top of the injector, undo and remove the screw closed to the coolant tube insertion hole. See [Figure 559](#).

**Figure 559.** Coolant Tube Insertion Hole

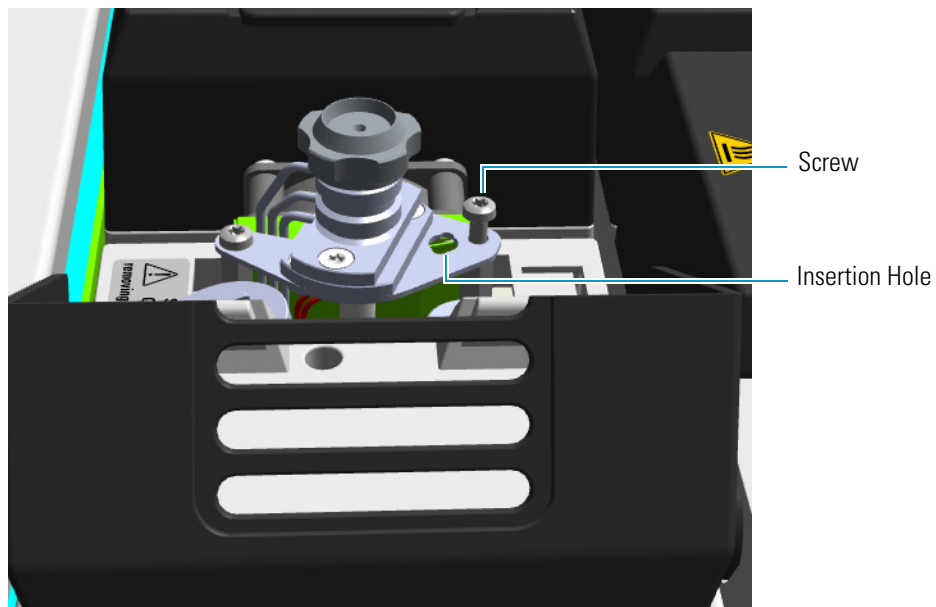

- c. Insert the coolant tube into the insertion hole, and fix it using the screw previously removed. See [Figure 560](#).

**Figure 560.** Installation of the Coolant Tubes for Carbon Dioxide and Liquid Nitrogen

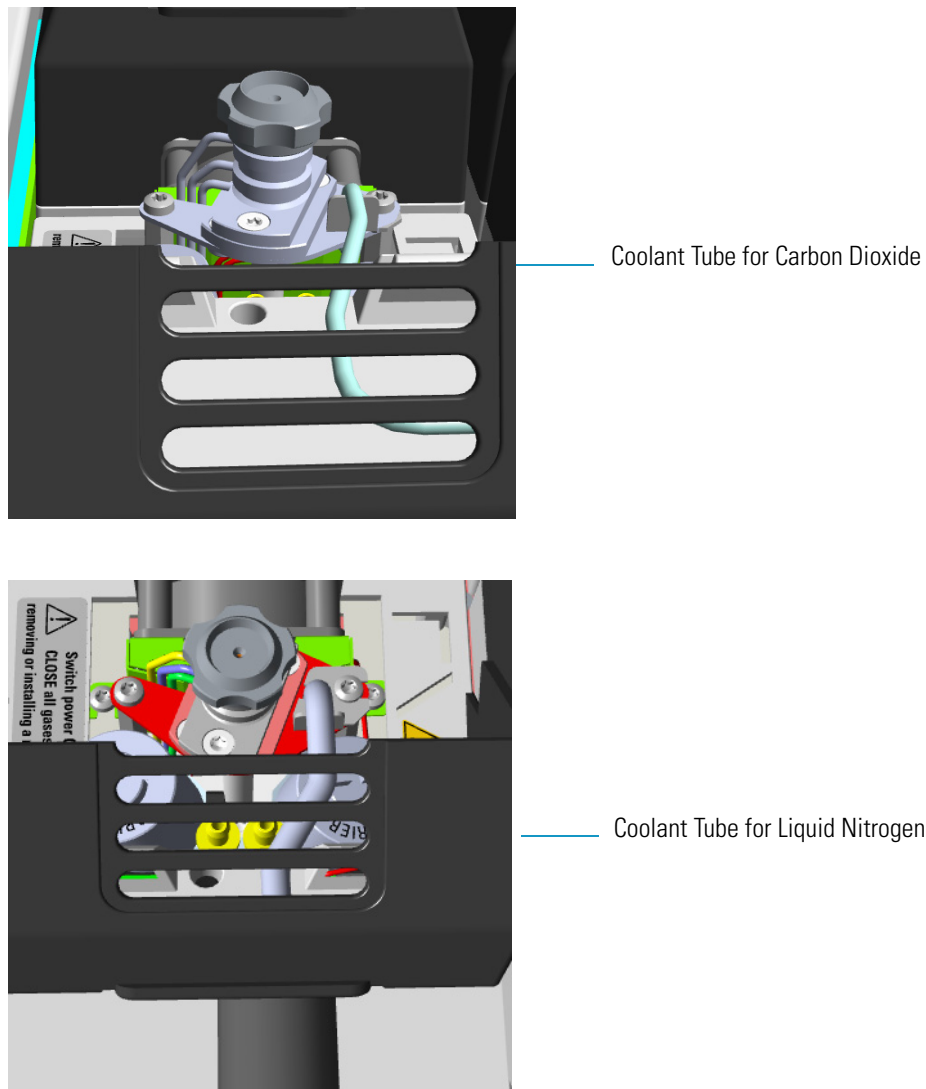


- d. Guide the coolant tube along the GC top cover up to reach the solenoid valve assemble on the back of the GC. Bend the tube if necessary.
11. Connect the coolant tube to the solenoid valve assemble.
    - a. Connect the coolant tube to the solenoid valve assembly using the proper nut and ferrule. Use a 7/16-in. wrench for tightening the fittings. See [Figure 561](#).

## 10 Adding Systems


Adding the PTV and PTVBKF Cryo System

**Figure 561.** Coolant Tube Connection to the Solenoid Valve

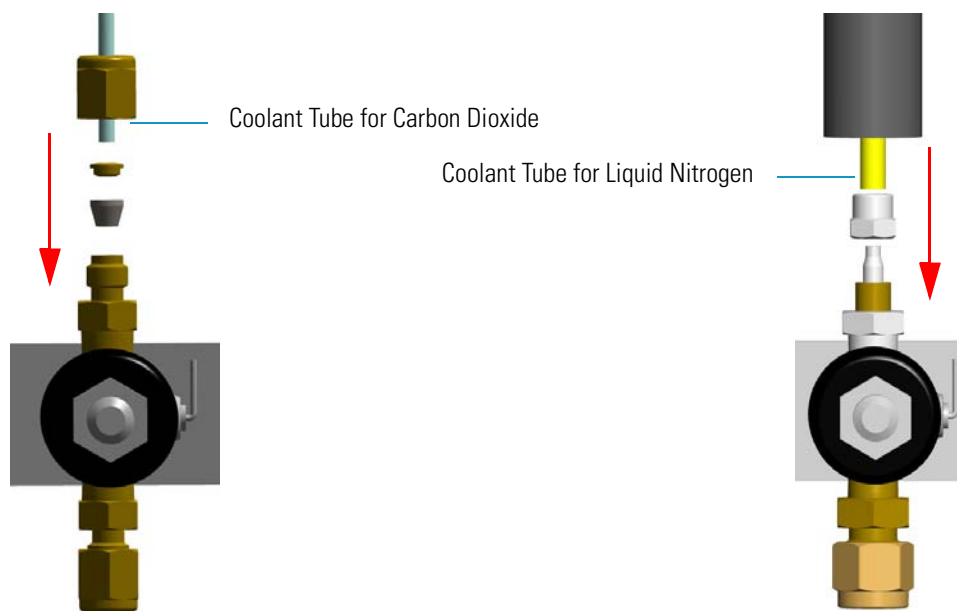


- b. Connect the proper end of the cryo supply tube to the solenoid valve using the proper nut and ferrule. Use a 7/16-in. wrench for tightening the fittings. See [Figure 562](#).

**Figure 562.** Cryo Supply Tube Connection to the Solenoid Valve

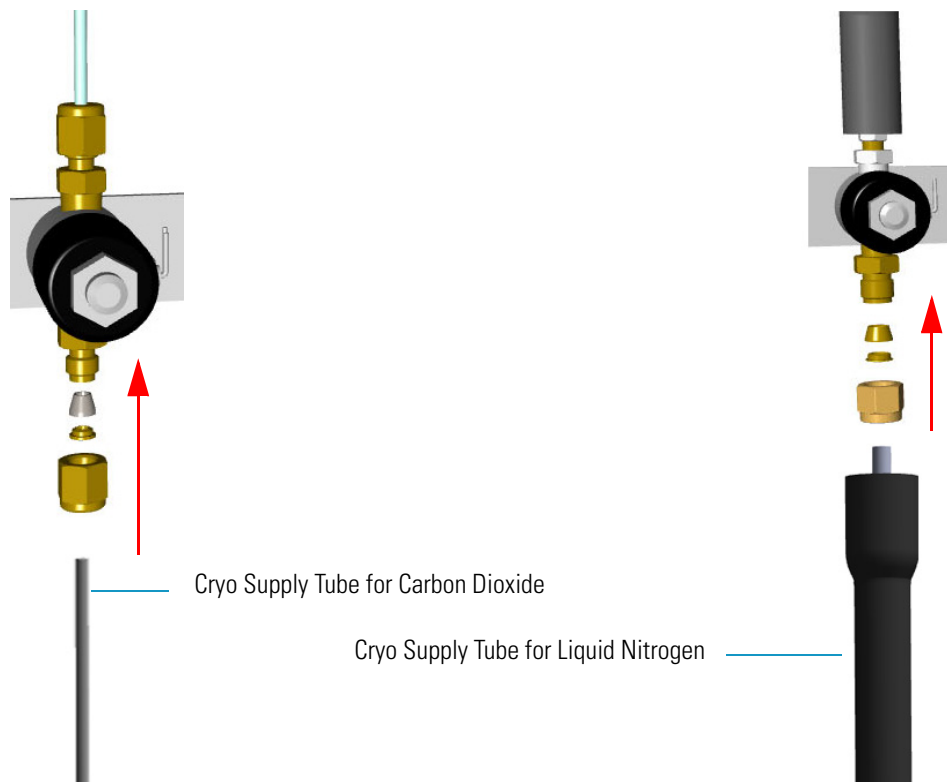

- c. Connect the other end of the cryo supply tube to the coolant container using the appropriate nuts and ferrules. Use a 7/16-in. wrench for tightening the fittings.
12. Connect the cryo valve to the Aux Temperature/Cryo Module.
  - a. Connect the cryo solenoid valve to the 2-pin connector marked **Cryo Valves - Front Inlet** or **Cryo Valves-Back Inlet** using the cable provided. See [Figure 563](#).

**Figure 563.** Cryo Valves: Oven

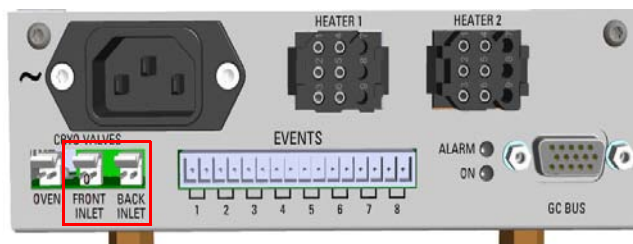

- b. Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.
  - c. Plug the power cable to the AC Input connector on the front of the module and to the wall outlet. The LED marked **On** lights after the GC is powered on.
- Note** For further details regarding the installation of the Aux Temperature/Cryo module, see the section [“Adding an Aux Temperature/Cryo Module”](#) on [page 443](#).
- d. Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.
  13. Mount the left side panel proceeding in reverse order which the left side panel was removed.
  14. Open the gas supplies.
  15. If other external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
  16. Power on the GC.
    - a. Plug the power cable to the AC Input connector on the back of the GC, and to the wall outlet.
    - b. Flip up the power switch (breaker), located at the back of the instrument, to the position I.
  17. Configure and enable the cryogenic system through the user interface of your GC, or through the CDS in use. Refer to the *TRACE 1300 and TRACE 1310 User Guide*.
    - a. **Configuration and enabling through the touch screen.**
      - i. In the main menu select the **Configuration** icon, the Configuration menu appears.

- ii. In the Configuration menu, select the **Front/Back PTV or PTVBKF** icon to open the relevant submenu.
    - iii. Set the **Cryogenic** parameters.
      - **Cryogenic Type** — Enable or disable the cryogenic system when it is installed and configured with Carbon Dioxide or Liquid Nitrogen as a coolant. Select between LN<sub>2</sub>, CO<sub>2</sub>, none.
      - **Cryo timeout** — Enter the time at which the cryo system will be disabled. Enter a value from 0–30 min.
      - **Cryo Threshold** — Specify the temperature at which the cryo system begins to supply the coolant. Enter a value from 40-200 °C.
      - **Cryo Cool at** — Specify the temperature at which the cryogenic system begins to supply the coolant. Enter a value from 40-200 °C.
    - iv. Return to main menu.
  - b. **Configuration and enabling through the Chromatography Data System.**
    - i. In the **Configuration** window select the **Auxiliary** tab.
    - ii. Select the **Auxiliary control module** check box to enable the setting for the auxiliary control of module and option installed on your GC.
    - iii. Select the **Front inlet cryogenics** or/and **Back inlet cryogenics** check box.
    - iv. Select the **Cryo type** used by your cryogenic option. Choose one: Liquid Nitrogen or Carbon Dioxide.
    - v. Open the **PTV/PTVBKF** page. Select the **Cryogenics enable** check box to enable the cryogenic system.
    - vi. In **Cool during** combo box, select when you want the cooling to be done. Choose between Prep-Run or Post -Run.
    - vii. In **Cryo threshold** text box specify the temperature at which the cryo system begins to supply the coolant. Enter a value in the range 40-200 °C.
    - viii. In **Cryo timeout** text box enter the time at which the cryo system will be disabled. Enter a value in the range 0–30 min.
18. If present, move the autosampler towards the module to restore the original alignment.
19. Set the normal detector, injector, and GC working conditions.

## Adding an Auxiliary Gas System

This section provides instructions for updating your TRACE 1300/TRACE 1310 with the Auxiliary Gas System.




**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.


The Auxiliary Gas Module is shipped with a protecting plate screwed on the manifold. Remove this plate before installing the module. See [Installing and Connecting the Auxiliary Gas Module](#) on page 534 for details.

The auxiliary gas system comprises the auxiliary gas module and the auxiliary gas interface. See [Figure 564](#).

**Figure 564.** Auxiliary Gas Module and Interface

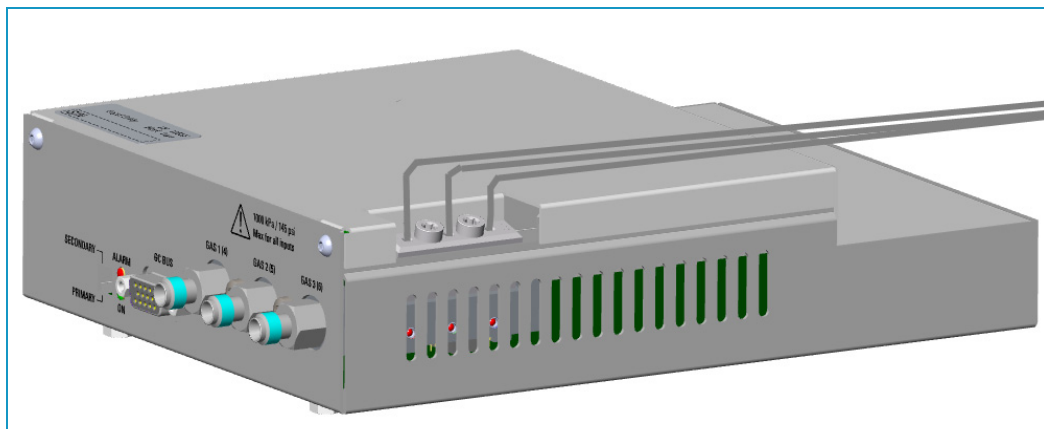

Auxiliary Gas Module with Gas Tubing Block

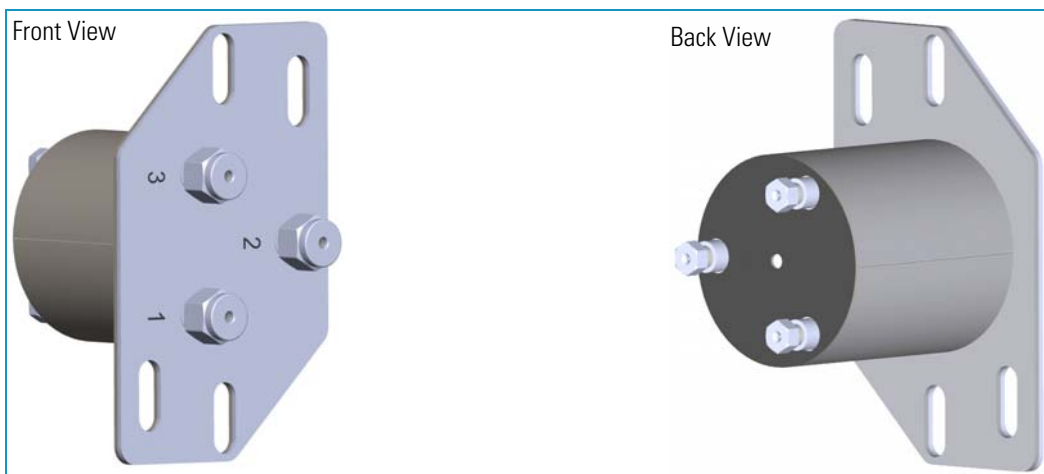


Auxiliary Gas Interface

## Auxiliary Gas Module Overview

The **auxiliary gas module** includes the following connections. See [Figure 565](#)

**Figure 565.** Auxiliary Gas Module Connections

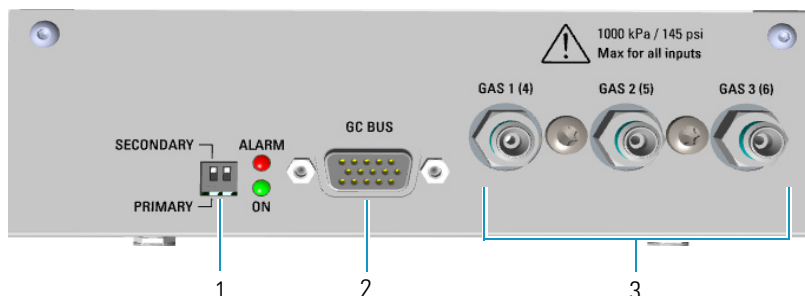


1. Switch marked **Primary/Secondary** used to set two Auxiliary Gas modules simultaneously present, one as Primary, and the other as Secondary. The primary module controls the aux pressures from 1 to 3, while the secondary module controls the aux pressures from 4 to 6. See also the point 3.
2. 15-pin female connectors marked **Bus** for the communication with the GC.
3. Three inlet ports marked **Gas 1 (4)**, **Gas 2 (5)**, and **Gas 3 (6)** for the connection up to three auxiliary carrier gases. If two modules are present, up to six auxiliary carrier gases can be connected. See also the point 1.

## Auxiliary Gas Interface Overview


The **auxiliary gas interface** is installed and fixed on the left or right wall of the GC oven, through the ducts provided, for the coupling with a mass spectrometer. See [Figure 566](#).

**Note** The Auxiliary Gas Interface is fixed on the exterior wall of the GC oven through the slots provided on the collar duct. Align the slots to the fixing holes accordingly.


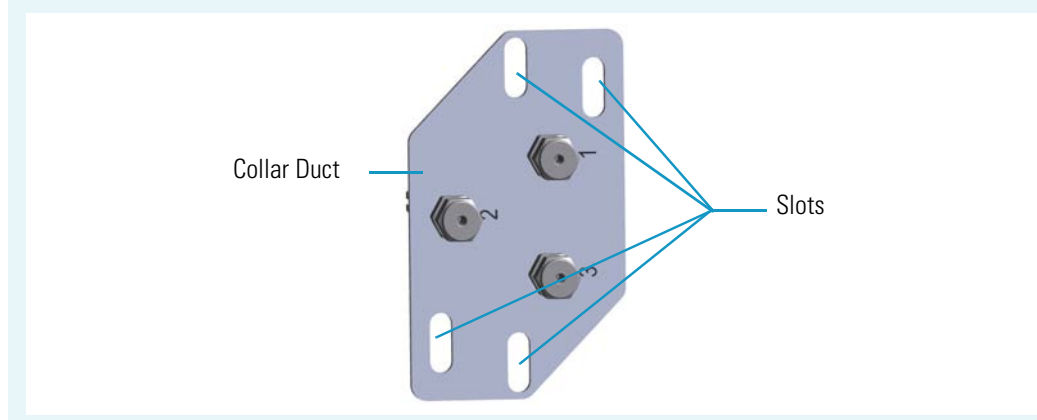





**Figure 566.** Auxiliary Gas Interface Installed in the Oven

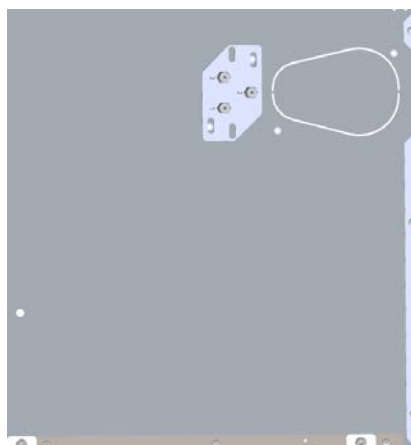


Auxiliary Gas Interface  
installed on the Left

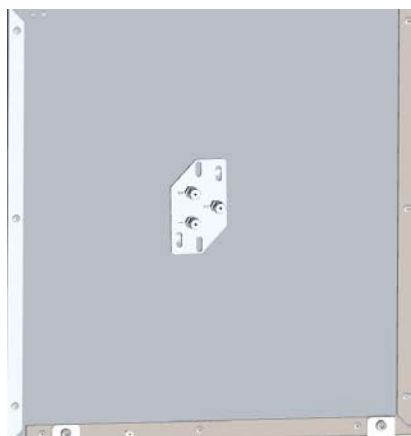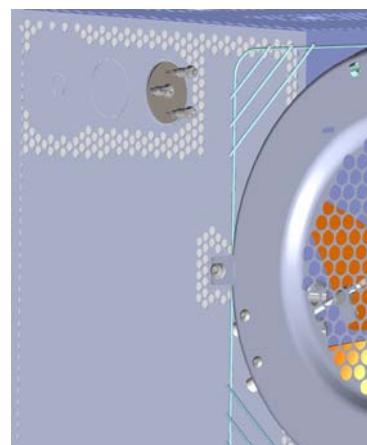

Auxiliary Gas Interface  
installed on the Right

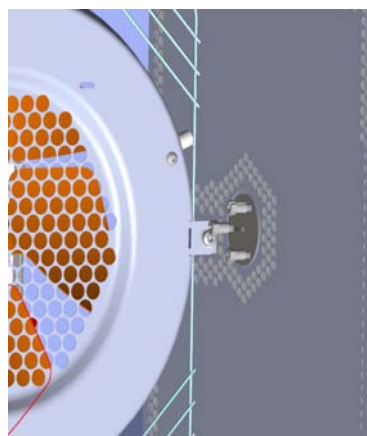


**Note** If your GC is equipped with the oven for the coupling with a Thermo Scientific high resolution mass spectrometers (HRMS), the Auxiliary Gas Interface is installed through the ducts provided on the left and right walls of the oven as well as the GC equipped with the standard oven. See [“Installing the Auxiliary Gas Interface on the Oven for HRMS”](#) on [page 532](#).

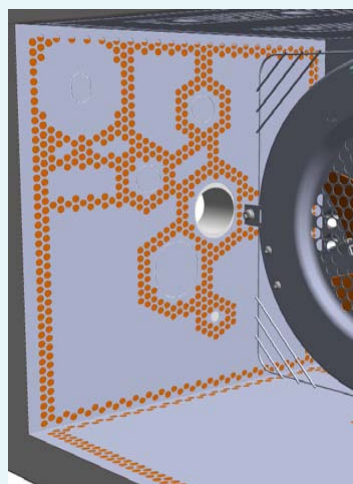


Left Wall

Right Wall

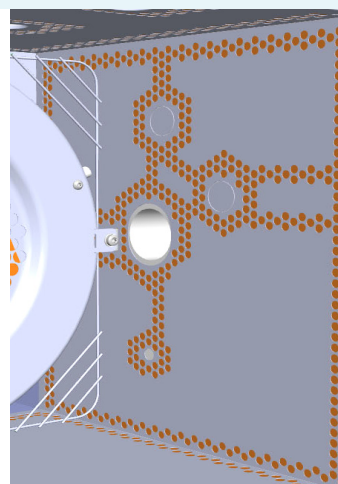

The module should be installed into a free external module housing provided on the back of the GC. See [Figure 567](#).

**Figure 567.** External Modules Housing

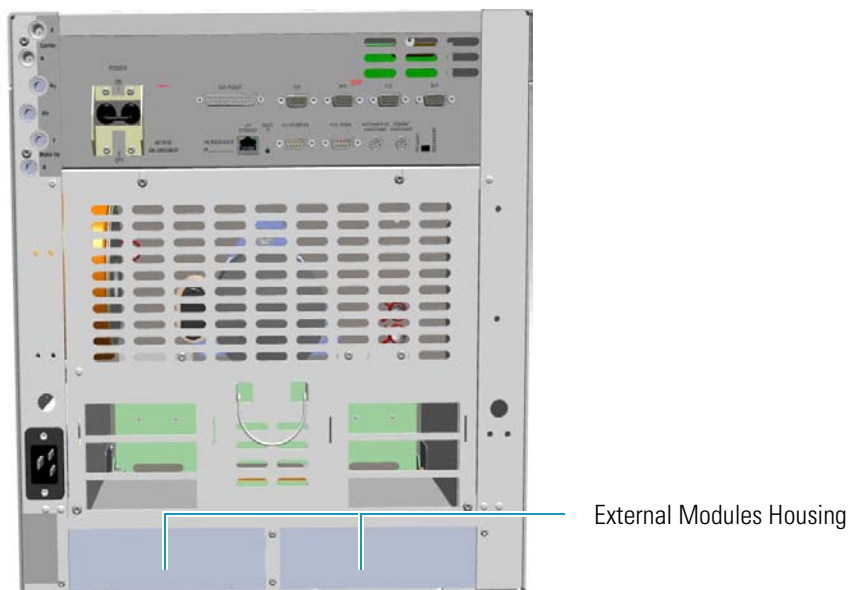

## Preliminary Operations

Before starting, the following preliminary operation must be carried out.

1. Put the GC in standby condition.
2. Cool the oven, injectors, and detectors to room temperature.

**Note** By pressing the **Maintenance** button, the GC cool down is automatically carried out.

3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector on the back of the GC, and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.

## Getting Started

To install and connect the Auxiliary Gas Interface and the Auxiliary Gas Module, see the following sections:


- “Installing the Auxiliary Gas Interface on the Left Wall of the Oven” on page 525
- “Installing the Auxiliary Gas Interface on the Right Wall of the Oven” on page 528
- “Installing the Auxiliary Gas Interface on the Oven for HRMS” on page 532
- “Installing and Connecting the Auxiliary Gas Module” on page 534

## Installing the Auxiliary Gas Interface on the Left Wall of the Oven

### ❖ To install the auxiliary gas interface on the left wall of the oven


1. Make sure that the preliminary operations have been carried out. See “[Preliminary Operations](#)” on page 524.
2. Remove the left side panel.
  - a. Open the front door of the GC.
  - b. Using a T20 Torxhead screwdriver, unscrew the left side panel screw from the interior front panel. See [Figure 568](#). Save the screw because it will be reused later.

**Figure 568.** Left Side Panel Fixing Screw

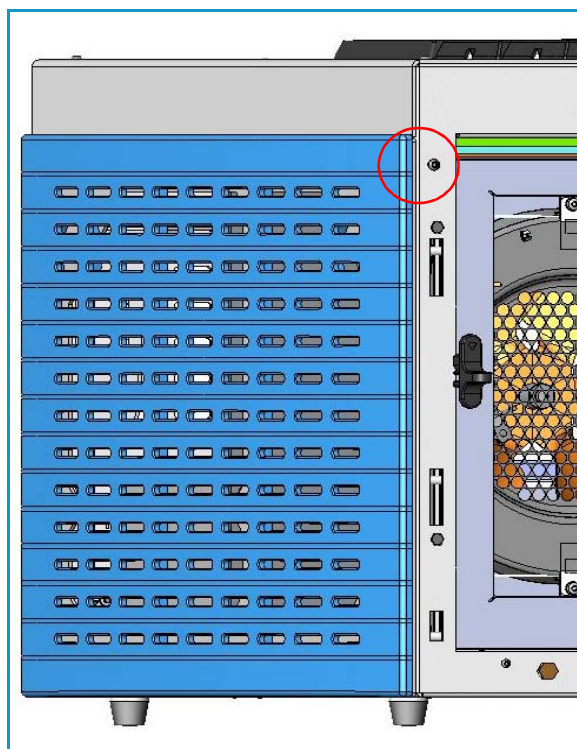


- c. Slide the panel towards the back of the instrument up to the stop.
- d. Remove the panel pulling it outward being aware that the ground wire is attached to the panel. See [Figure 569](#).
- e. Unplug the ground wire from the panel.

**Figure 569.** Left Side Panel Removal

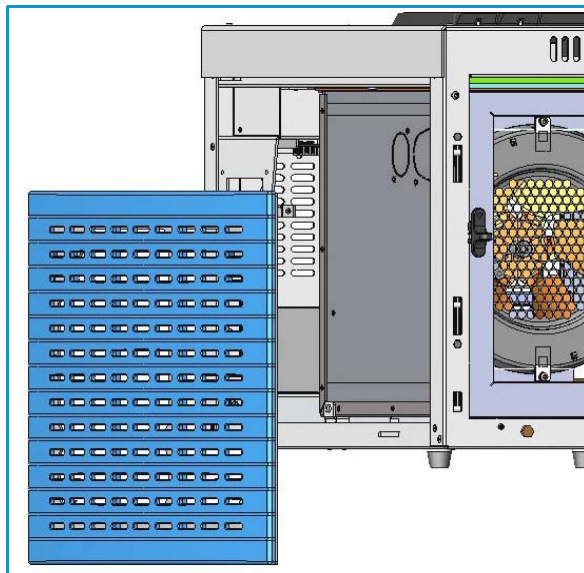


3. Prepare the duct for the installation of the auxiliary gas interface.
  - a. Remove the partial cut shaped plate from the exterior wall of the oven box for accessing the insulating material. See [Figure 570](#).

**Figure 570.** Perform the Duct for the Auxiliary Gas Interface (1)

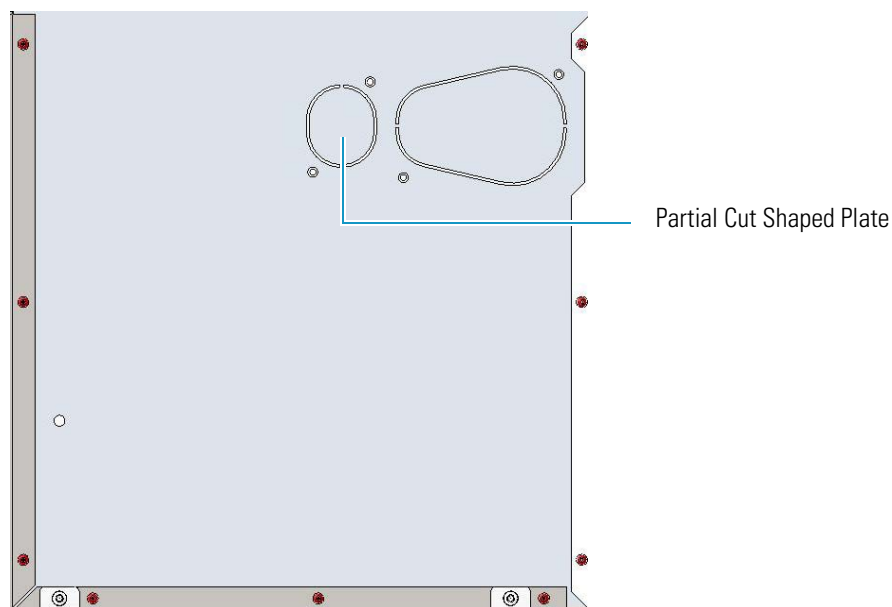


- b. Using a knife or similar tool, gently cut the insulating material following the border. See [Figure 571](#).

**Figure 571.** Perform the Duct for the Auxiliary Gas Interface (2)

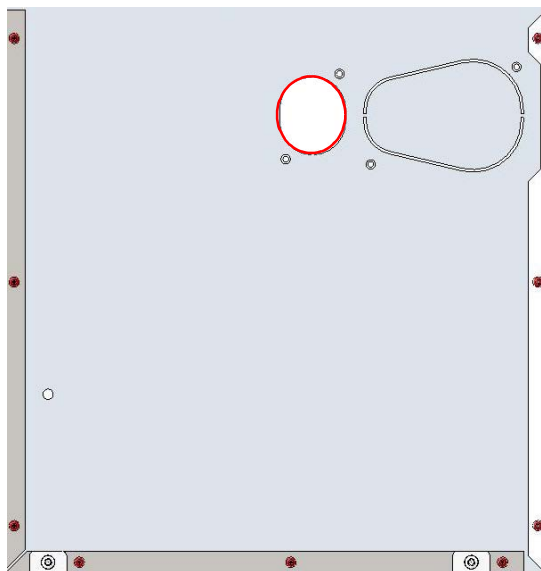


- c. Save the removed insulating material in a safe place because it can be reused.
- d. On the left side wall in the interior of the oven, remove the partial cut plate from the corresponding duct. See [Figure 572](#).

**Figure 572.** Perform the Duct for the Auxiliary Gas Interface (3)

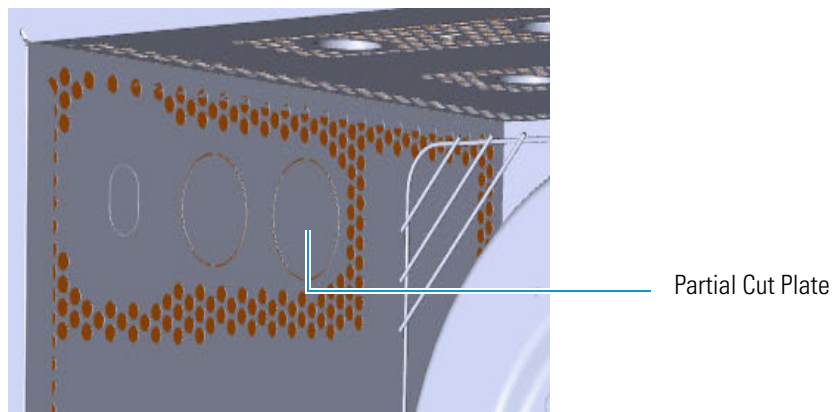

- e. Insert the auxiliary gas interface into the duct. Fix the interface on the exterior wall of the oven box using the fixing screws provided. See [Figure 573](#).

**Figure 573.** Perform the Duct for the Auxiliary Gas Interface (4)

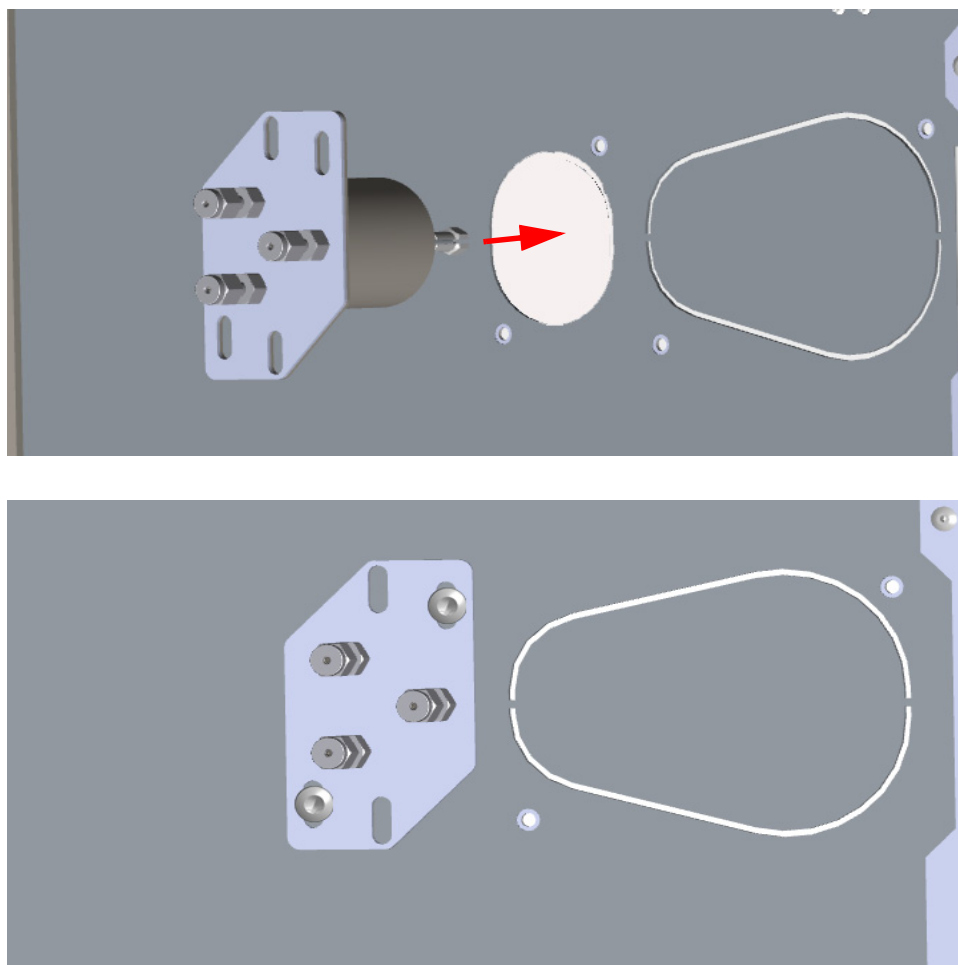


4. Jump to the section [“Installing and Connecting the Auxiliary Gas Module”](#) on [page 534](#).

## Installing the Auxiliary Gas Interface on the Right Wall of the Oven

### ❖ To install the auxiliary gas interface on the right wall of the oven


1. Make sure that the preliminary operations have been carried out. See [“Preliminary Operations”](#) on [page 524](#).
2. Remove the right side panel.
  - a. Open the front door of the GC.
  - b. Using a T20 Torxhead screwdriver, unscrew the left side panel screw from the interior front panel. See [Figure 574](#). Save the screw because it will be reused later.

**Figure 574.** Right Side Panel Fixing Screw

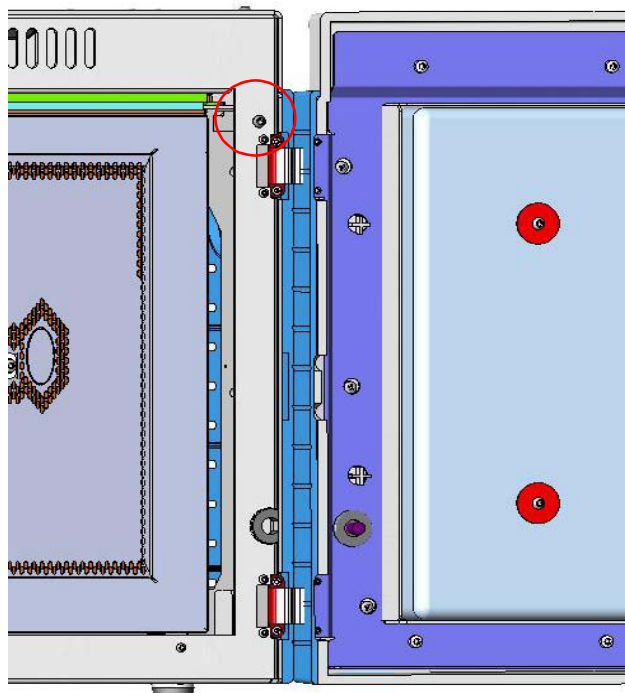

- c. Slide the panel towards the back of the instrument up to the stop.
- d. Remove the panel pulling it outward. Be aware that the ground wire is attached to the panel. See [Figure 575](#).
- e. Unplug the ground wire from the panel.

**Figure 575.** Right Side Panel Removal

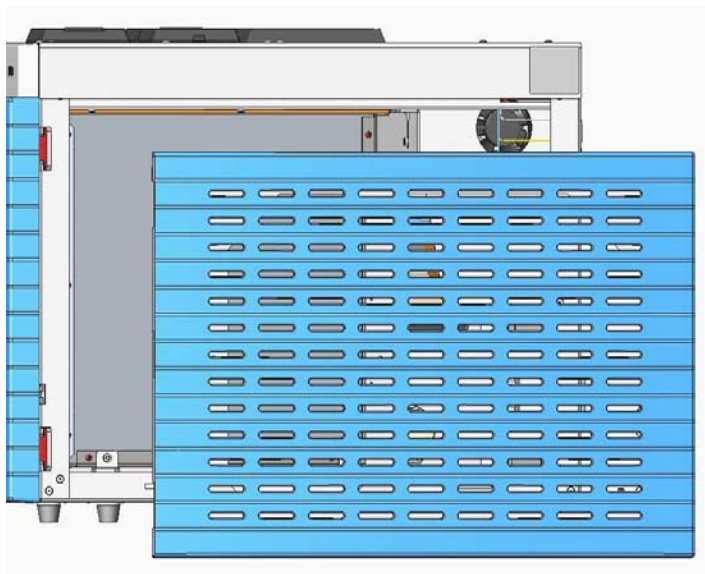


3. Prepare the duct for the installation of the auxiliary gas interface.




- a. Remove the partial cut shaped plate from the exterior wall of the oven box for accessing the insulating material. See [Figure 576](#).

**Figure 576.** Perform the Duct for the Auxiliary Gas Interface (1)

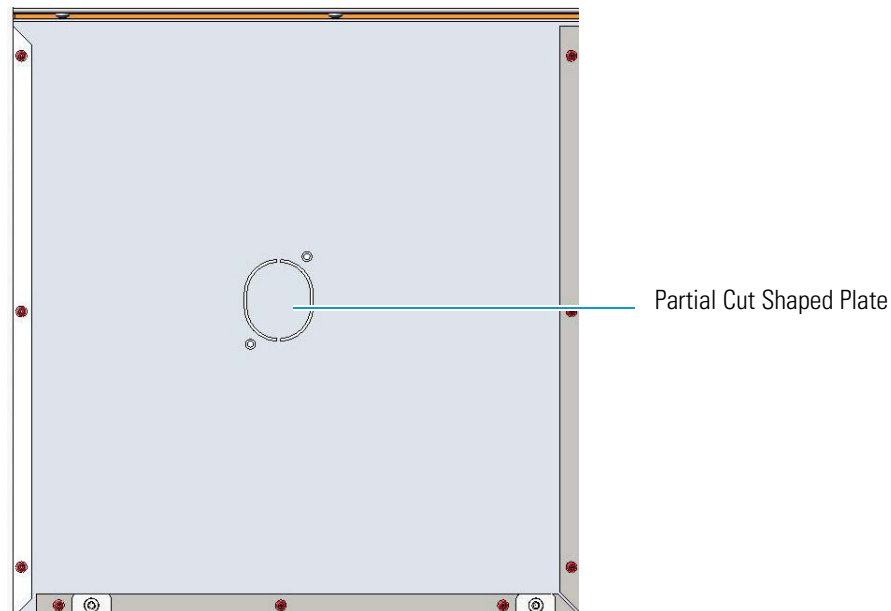

- b. Using a knife or similar tool, gently cut the insulating material following the border. See [Figure 577](#).


**Figure 577.** Perform the Duct for the Auxiliary Gas Interface (2)

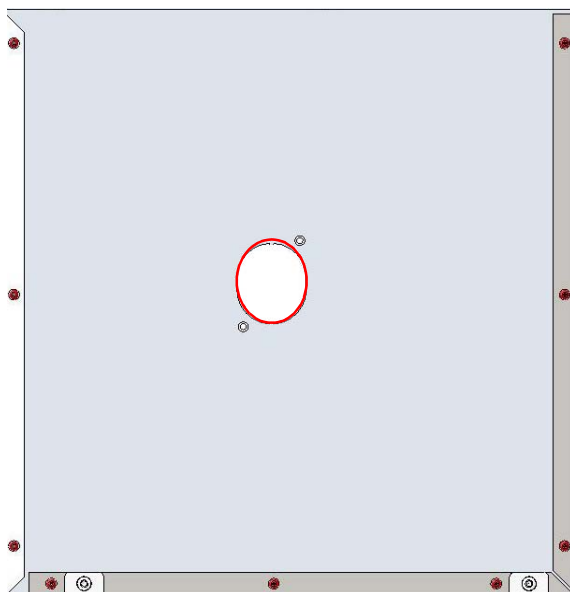

- c. Save the removed insulating material in a safe place because it can be reused.
- d. On the left side wall in the interior of the oven, remove the partial cut plate from the corresponding duct. See [Figure 578](#).




**Figure 578.** Perform the Duct for the Auxiliary Gas Interface (3)

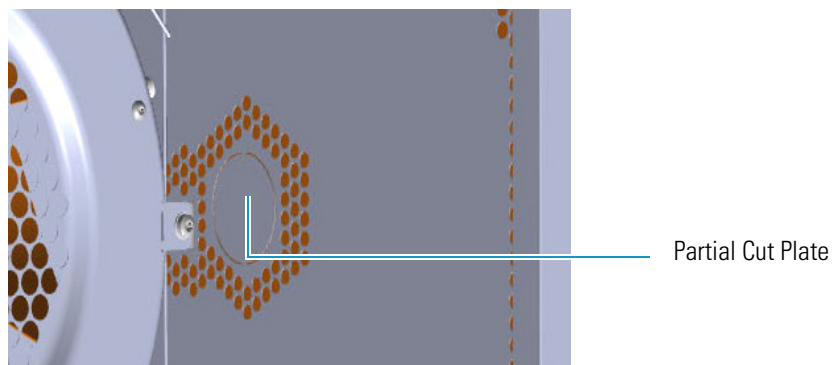

- e. Insert the auxiliary gas interface into the duct. Fix the interface on the exterior wall of the oven box using the fixing screws provided. See [Figure 579](#).

**Figure 579.** Perform the Duct for the Auxiliary Gas Interface (4)

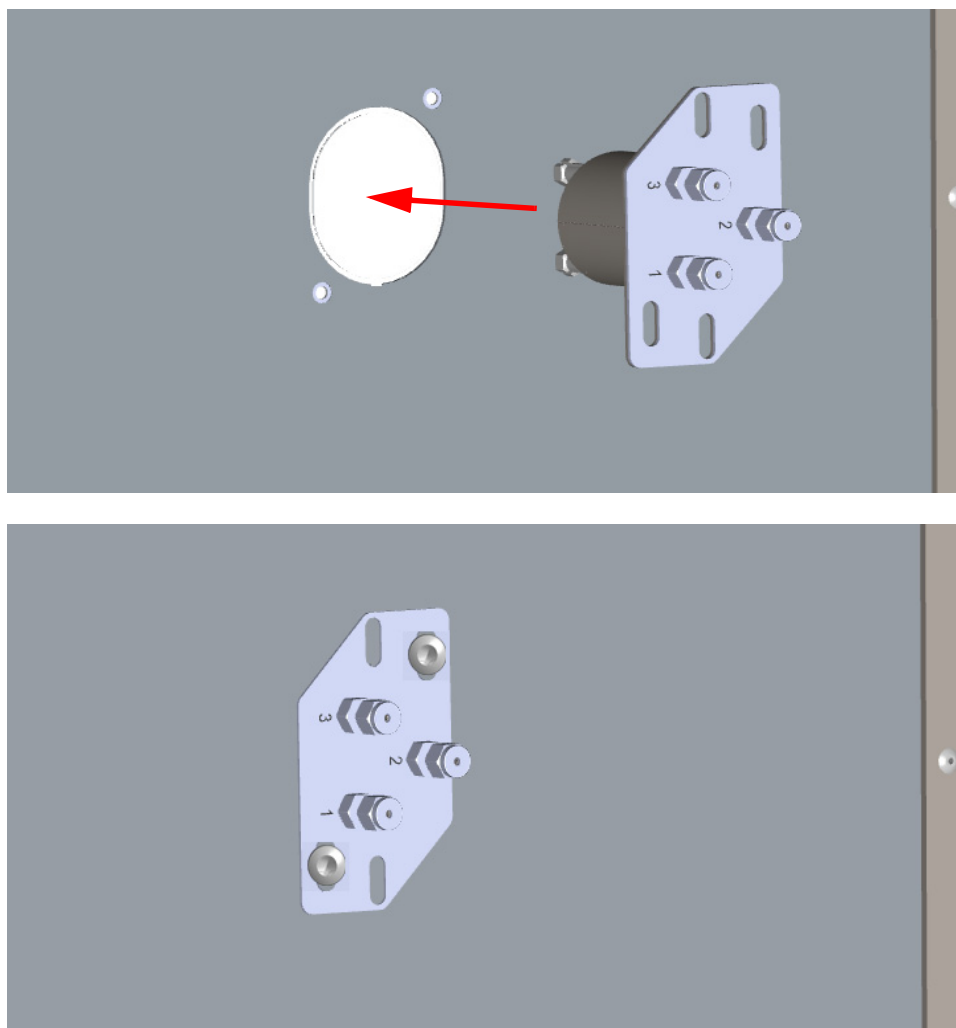


- 4. Jump to the section “[Installing and Connecting the Auxiliary Gas Module](#)” on [page 534](#).

## Installing the Auxiliary Gas Interface on the Oven for HRMS

### ❖ To make the duct into the oven for HRMS

1. Remove the left/right side panel.
  - a. Open the front door of the GC.
  - b. Use a T20 Torxhead screwdriver to loosen the left/right side panel screw from the interior front panel. Save the screw because it will be reused later.
  - c. Slide the panel towards the back of the instrument up to the stop.
  - d. Remove the panel pulling it outward being aware that the ground wire is attached to the panel.
2. Prepare the duct for the installation of the auxiliary gas interface.
  - a. On the left/right exterior wall of the oven box, remove the partial cut plate of the duct of interest.
  - b. Remove the pre-shaped plug of insulating material from the duct provided. See [Figure 580](#).


**Figure 580.** Left/Right Exterior Oven Wall Box View

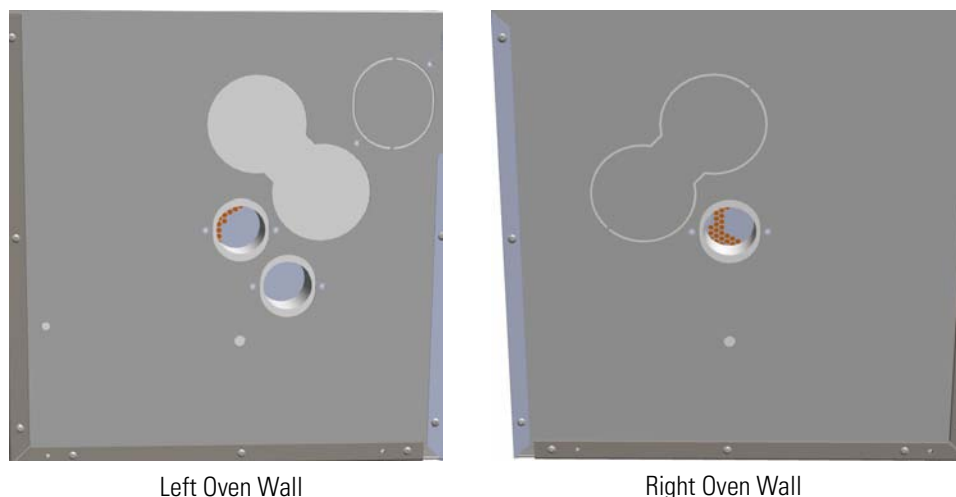

**IMPORTANT** Save the pre-shaped plug of insulating material in a safe place because it could be reused.

- c. On the interior of the oven box, remove the partial cut plate from the corresponding duct. See [Figure 581](#).

**Figure 581.** Left/Right Interior Oven Wall Box View

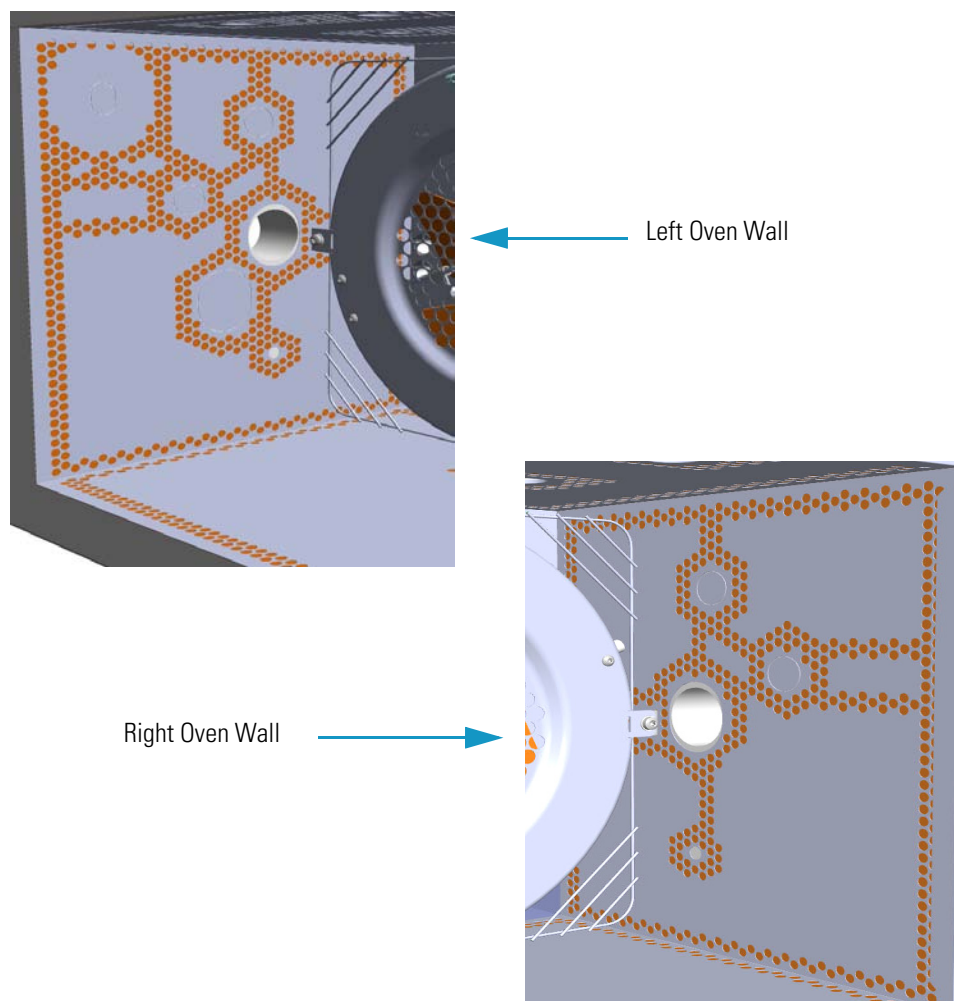


3. Insert the auxiliary gas interface into the duct. Fix the interface on the exterior wall of the oven box using the fixing screws provided.
4. Jump to the section [“Installing and Connecting the Auxiliary Gas Module”](#) on [page 534](#).

## Installing and Connecting the Auxiliary Gas Module

### ❖ To install and connect the Auxiliary Gas Module


1. Remove the cover of the external modules housing where installing the module. See [Figure 582](#).

**Figure 582.** Housing Cover Removal

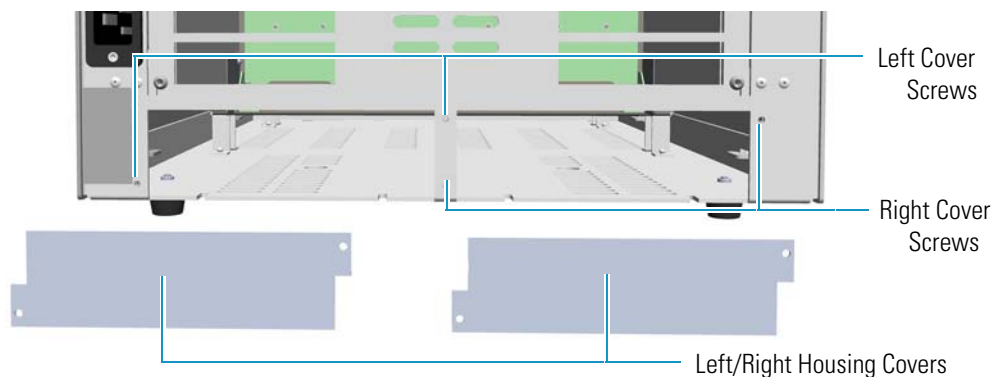

- a. Using a T20 Torxhead screwdriver, unscrew and remove the left or right housing cover screws.
  - b. Remove the cover from the housing.
2. Remove the manifold protecting plate
    - a. Using a T20 Torxhead screwdriver, unscrew the two fixing screws, and remove the protecting plate from the manifold. Save the protecting plate and the fixing screws.


**Figure 583.** Manifold Protecting Plate

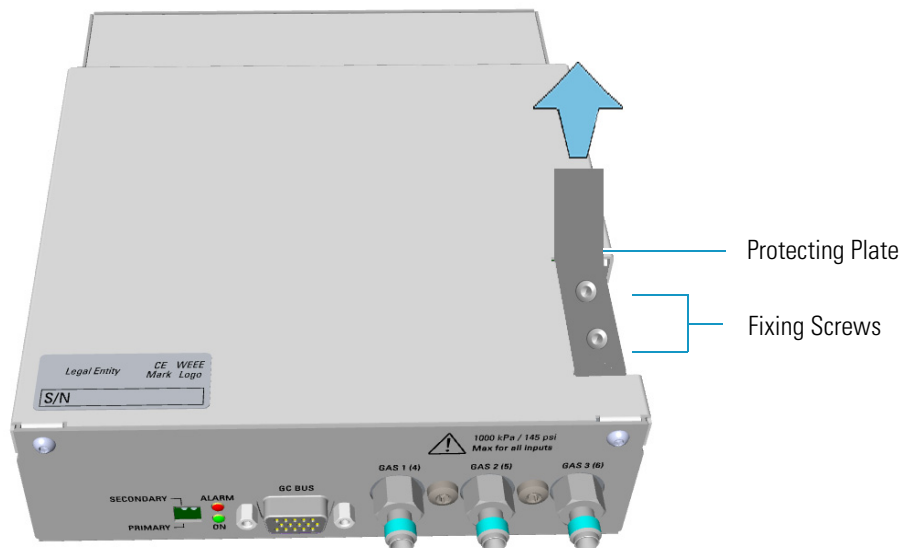

3. Install the module into the housing


- a. Loosen the two hexagonal screws under the module. See [Figure 584](#).

**Figure 584.** Module Installation (1)

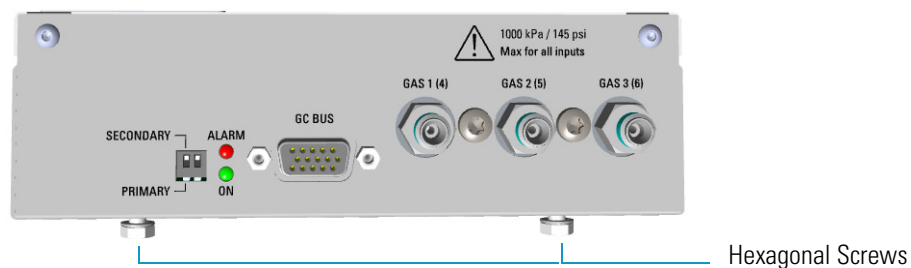

- b. Carefully place the module into the left or right housing.
- c. Push the module until the hexagonal screws couple with the slots on the floor of the GC. See [Figure 585](#).

**Figure 585.** Module Installation (2)

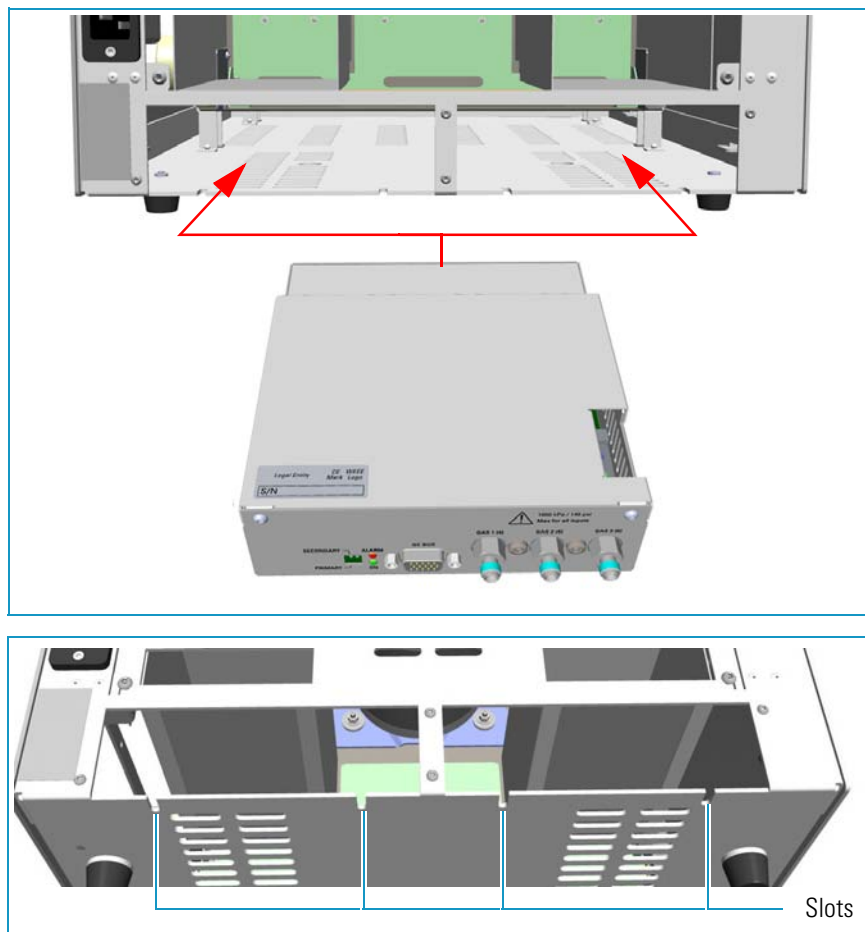


- d. Finger-tighten the hexagonal screws slightly, or use a 10-mm wrench.

**Note** Always keep the hexagonal screws in their place. This allows you an easier removal of the auxiliary module when necessary.

The result of the installation is shown in [Figure 586](#).


**Figure 586.** Auxiliary Gas Module Installed into the GC

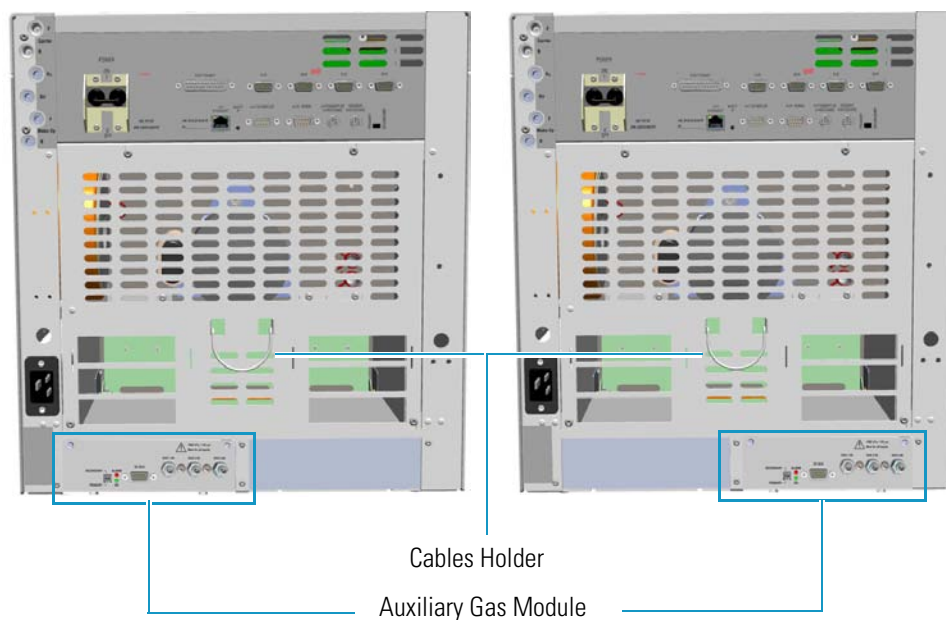

4. Connect the gas tubing block to the manifold.


[Figure 587](#) shows the gas tubing block and the manifold located into the auxiliary gas module.

**Figure 587.** Gas Tubing Block and Manifold

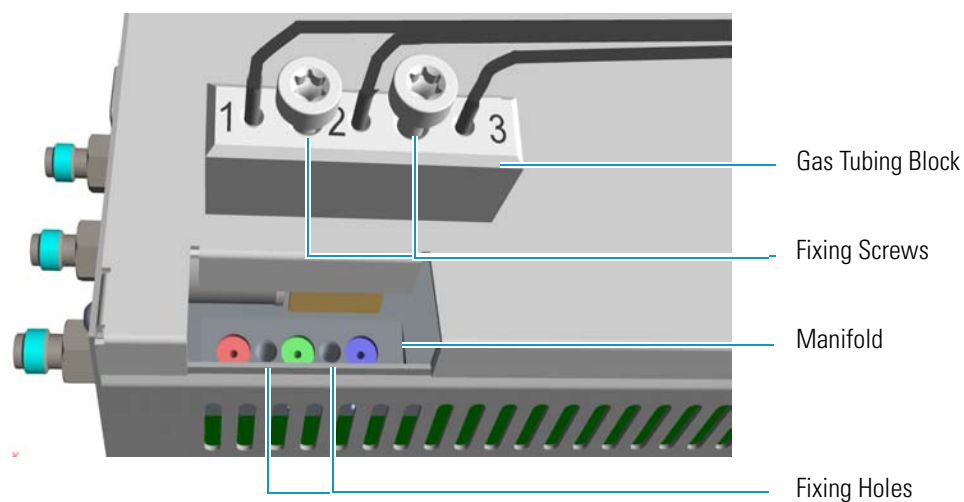


- a. Carefully guide the gas tubing block on the manifold located into the auxiliary gas module. See [Figure 588](#) if the module is installed on the right, or [Figure 589](#) if the module is installed on the left.

**Figure 588.** Installation on the Right Side

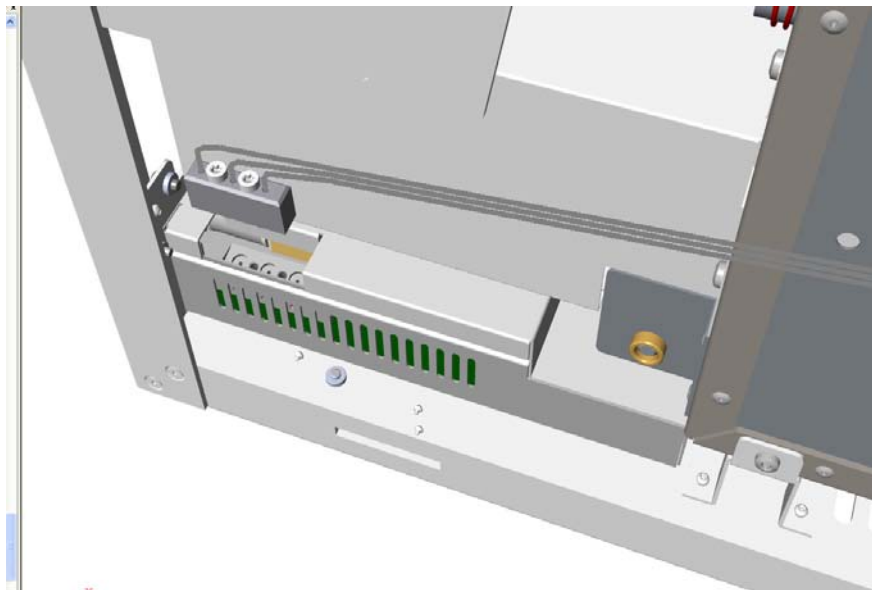


**Figure 589.** Installation on the Left Side

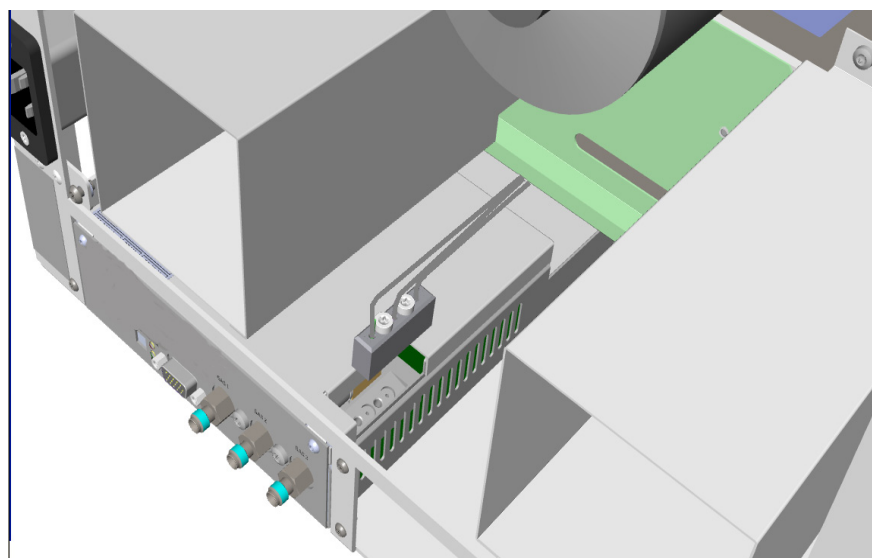

- b. Align the fixing screws of the gas tubing block with the corresponding holes on the manifold.
  - c. Use the T20 Torxhead screwdriver to tighten the two fixing screws without overtightening.
5. Connect the gas tubes of interest to the Auxiliary Gas Interface.
    - a. Guide the three gas tubes up to reach the Auxiliary Gas interface.


**Note** The length of the tubes allows them to reach the Auxiliary Gas Interface whether they are installed on the same side or on the opposite side of the Auxiliary Gas Module.

**ATTENTION** The three gas tubes, coming from the gas tubing block, are numbered 1, 2, and 3 respectively. Pay attention to the correct order when you connect each tube to the corresponding inlet on the auxiliary gas interface. The end of each tube is provided with a label indicating the type of gas.

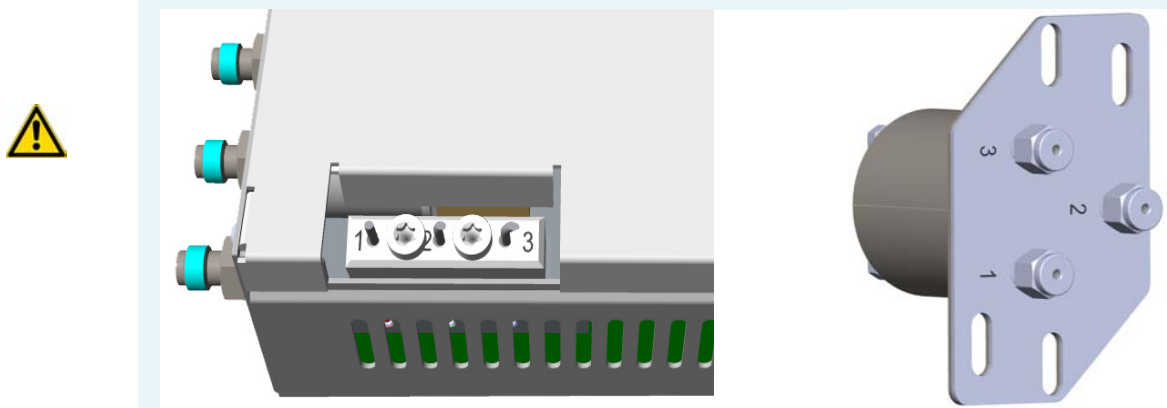

- b. Bend the gas tube until its end reaches the corresponding numbered inlet port of the auxiliary interface.
- c. Connect the gas tube to the corresponding numbered inlet port using the appropriate nut and ferrules. Use a 7/16-in. wrench to tighten the fittings. See [Figure 590](#).

**Figure 590.** Gas Tubes Connection to the Auxiliary Gas Interface - External View(

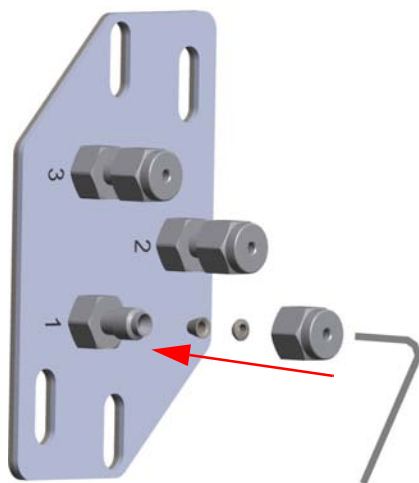


- d. Repeat [step b](#) and [step c](#) until all the gas tubes of interest are connected to the auxiliary gas interface.
- e. In the GC oven, carry out the connections of the components of interest to the corresponding inlet ports of the auxiliary gas interface using the appropriate nut and ferrules.



- f. Repeat [step e](#) until all the components are connected to the auxiliary gas interface. See [Figure 591](#).

**Figure 591.** Gas Tubes Connection to the Auxiliary Gas Interface - Internal View

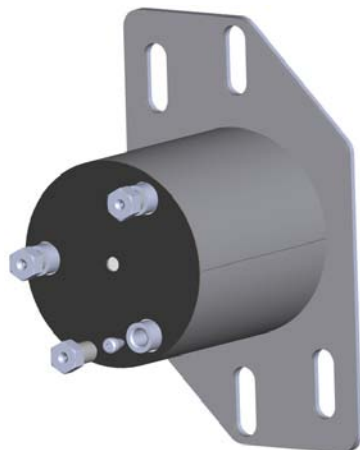


6. Connect the Supply Gas to the Auxiliary Gas Module.

**Note** Use the 1/8-in. Swagelok fittings provided on the gas inlet ports to connect the gas lines.

- a. Connect the gas line to the corresponding inlet port of interest using the appropriate nut and ferrules. Use a 7/16-in. wrench for tightening the fittings.
- b. Repeat [step a](#) until all the gas lines of interest are connected to the corresponding inlet port on the auxiliary gas module. See [Figure 592](#).

**Figure 592.** Gas Line Connection to the Auxiliary Gas Interface

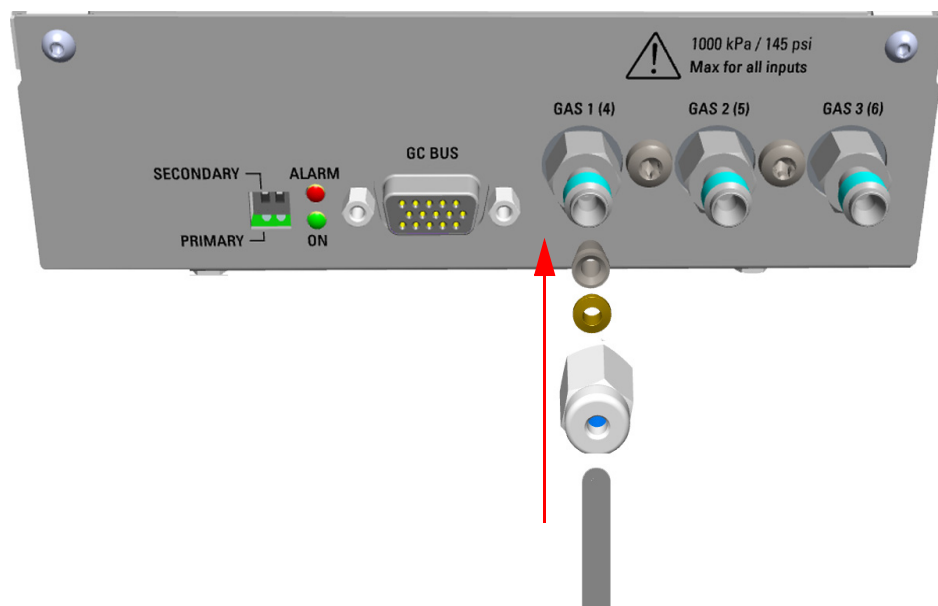




**IMPORTANT** The maximum nominal inlet pressure for all the inputs is 1000 kPa (145 psig). The working inlet pressure range is from 400 kPa (58 psig) to 1000 kPa (145 psig).

7. Connect the Auxiliary Gas Module electrically.
  - a. Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.

**IMPORTANT** If two Auxiliary Gas modules are simultaneously present, one must be set as Primary and the other as Secondary by mean of dedicated switch on the back panel. See the example below.

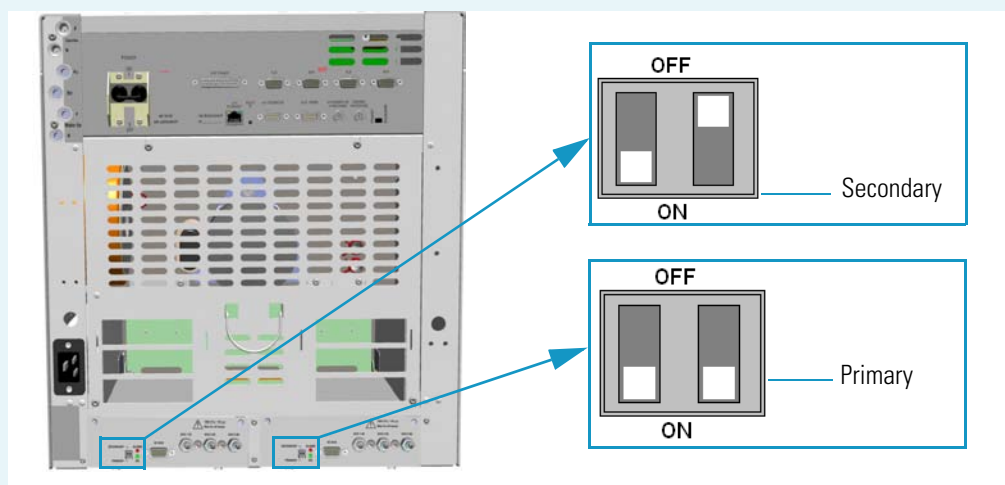

The primary module controls the aux pressures from 1 to 3, while the secondary module controls the aux pressures from 4 to 6.


8. Replace the left/right side panel.
  - a. Plug the ground wire previously removed into the left/right panel.
  - b. Place the left/right panel and attach the screw holding it in place.
9. Open the gas supplies.
10. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
11. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker), located at the back of the instrument, to the position I.

12. Configure and enable the Auxiliary Gas system through the user interface of your GC, or through the CDS in use. Refer to the *TRACE 1300 and TRACE 1310 User Guide*.
  - a. **Configuration and enabling through the touch screen.**
    - i. In the main menu select the **Instrument control** icon. The Instrument Control menu appears.
    - ii. In the Instrument Control menu, select the **Auxiliary** icon to open the relevant submenu.
    - iii. Set the **Aux Gas Pressure** values as required, then return to main menu.
  - b. **Configuration and enabling through the Chromatography Data System.**
    - i. In the **Configuration** window select the **Auxiliary** tab.
    - ii. Select the check box **Auxiliary Carrier Module 1/2** according to the auxiliary carrier module installed on your GC.
    - iii. Select the **Auxiliary Pressure** check box to enable up to six auxiliary Pressures, and the adjacent field.
    - iv. According the inlet ports connected to the Auxiliary Gas Interface, select the corresponding check box and set the pressure in the adjacent field.
13. Set the normal detector, injector, and GC working conditions.

## Adding the Hydrogen Sensor

This section provides instructions for updating your TRACE 1300/TRACE 1310 with the hydrogen sensor. See [Figure 593](#).

**Figure 593.** Hydrogen Sensor

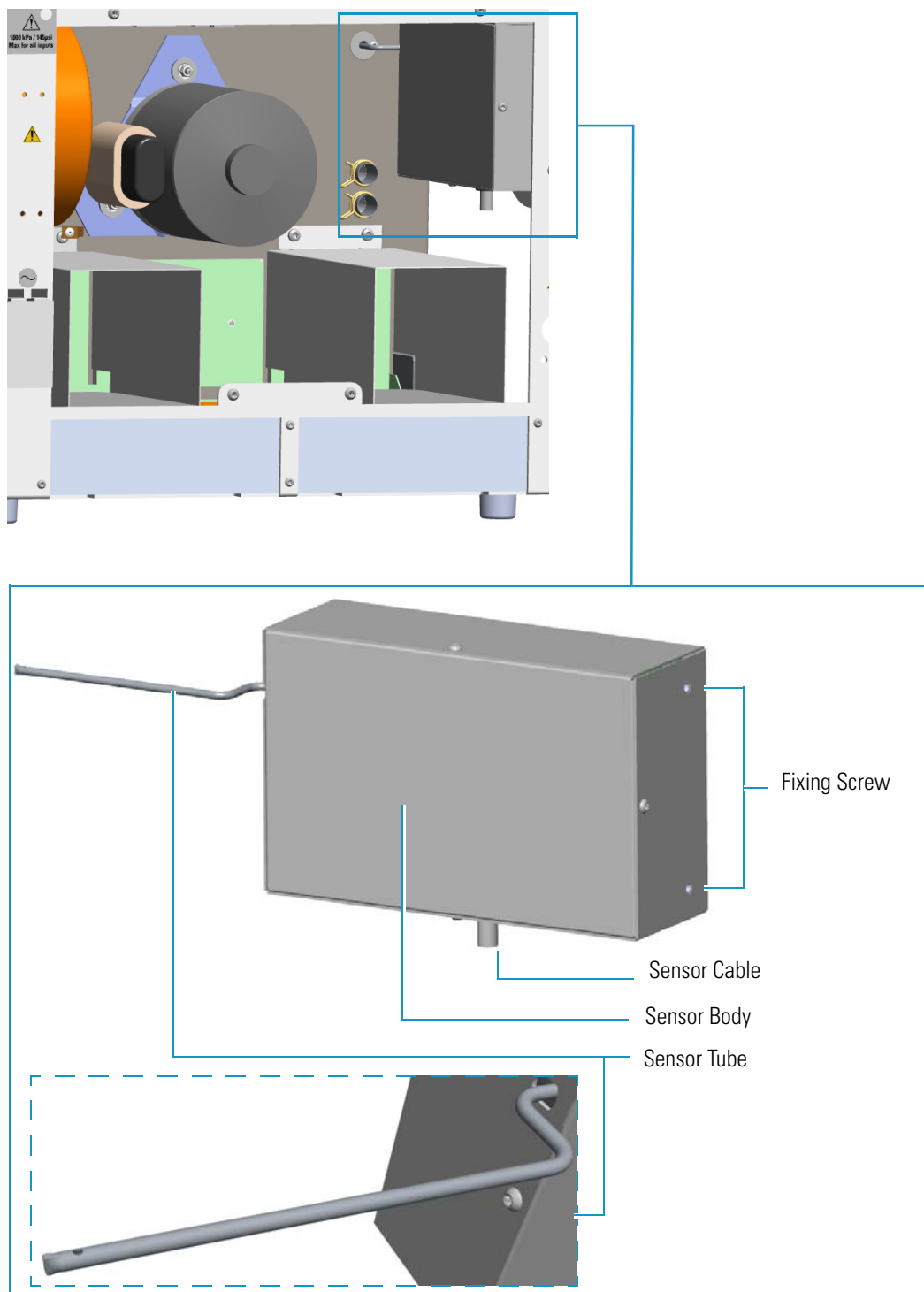

❖ **To add the hydrogen sensor**




**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.

**IMPORTANT** The hydrogen sensor requires a GC firmware version 1.03 or later.


1. Put the GC in standby condition.
  2. Cool the oven, injectors and detectors to room temperature.
  3. Close the gas supplies.
  4. Power off the GC.
    - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
    - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
  5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.
  6. Remove the left side panel.
    - a. Open the front door of the GC.
    - b. Using a T20 Torxhead screwdriver, unscrew the left side panel screw from the interior front panel. Save the screw because it will be reused later.
    - c. Slide the panel towards the back of the instrument up to the stop.
    - d. Remove the panel by pulling it outward. Be aware that the ground wire is attached to the panel.
    - e. Unplug the ground wire from the panel.
  7. Remove the back cover.
    - a. Use a T20 Torxhead screwdriver to remove the four screws that secure the back cover to the GC.
    - b. Lift the cover off using the cover handle. Be aware that the ground wire is attached to the back panel.
- Note** Pay attention to the positioning of the ground wire plug, so it can be reconnected in the same way it was removed.
8. Perforate the duct for the installation of the sensor into the oven.
    - a. Looking the interior of the GC from the back side, locate the duck provided on the back wall of the oven. See [Figure 594](#).

**Figure 594.** Duck for Hydrogen Sensor

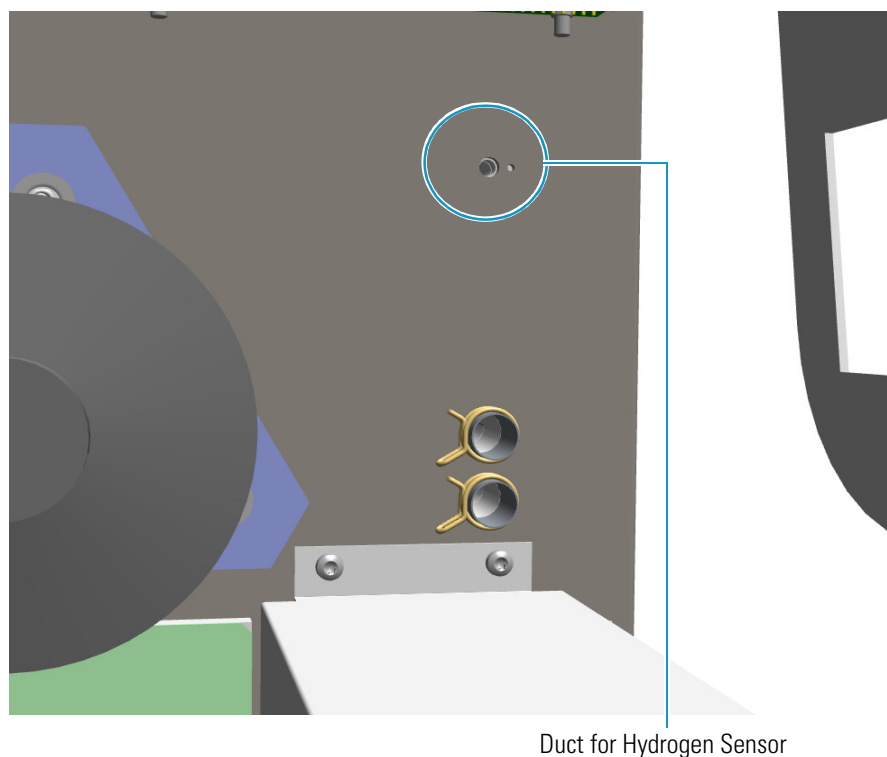


- b. Using a tool (for example a punch), perforate the insulating material until the tool protrudes into the oven. Make sure the duct is free of insulating material.
- c. Insert the sensor tube adapter into the duct and fix the adapter using the screw provided. See [Figure 595](#).

**Figure 595.** Sensor Tube Adapter Installation (1)

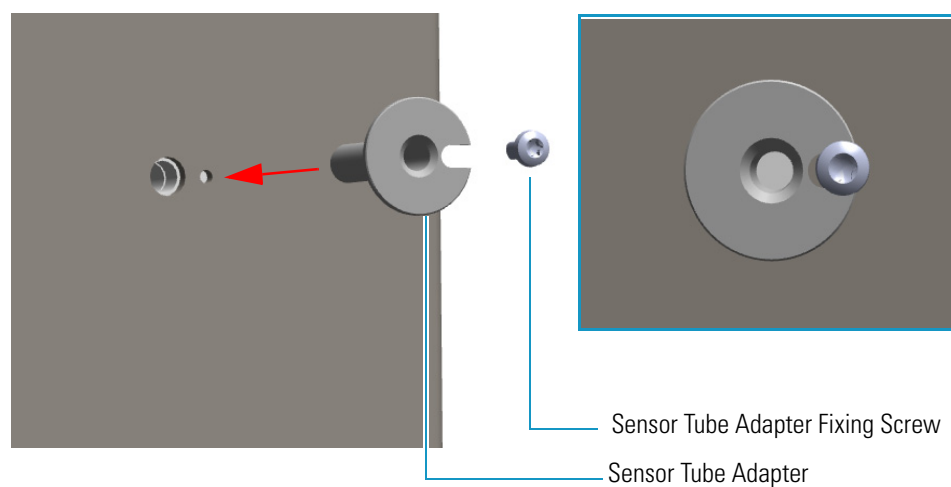

9. Open the front door.
  - a. Look into oven the duct for the hydrogen sensor previously done. See [Figure 596](#).

**Figure 596.** View of the Sensor Tube Duct Into the Oven

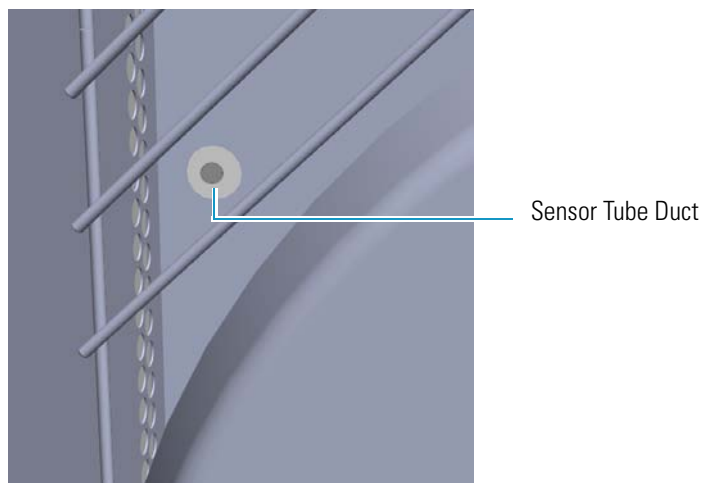

- b. Check the duct is free of insulating material; if not, remove it.




**IMPORTANT** The duct must be free of the insulating material. If not, it could obstruct the tube sensor.

10. Install the hydrogen sensor.


- a. Move and guide the sensor tube into the duct for the hydrogen sensor. See [Figure 597](#).

**Figure 597.** Hydrogen Sensor Installation (1)

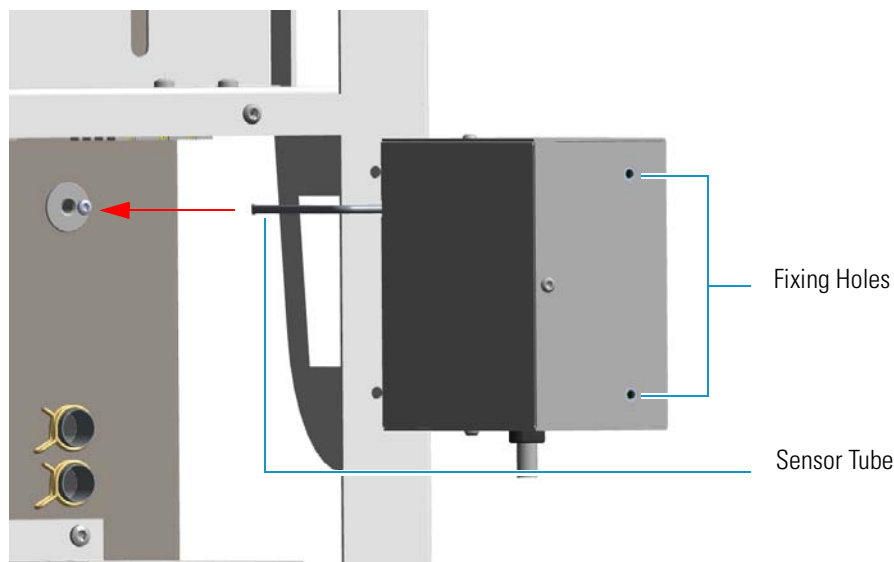


- b. Place the hydrogen sensor into the back of the GC aligning the fixing holes to the corresponding holes on the GC chassis. See [Figure 598](#).

**Figure 598.** Hydrogen Sensor Installation (2)

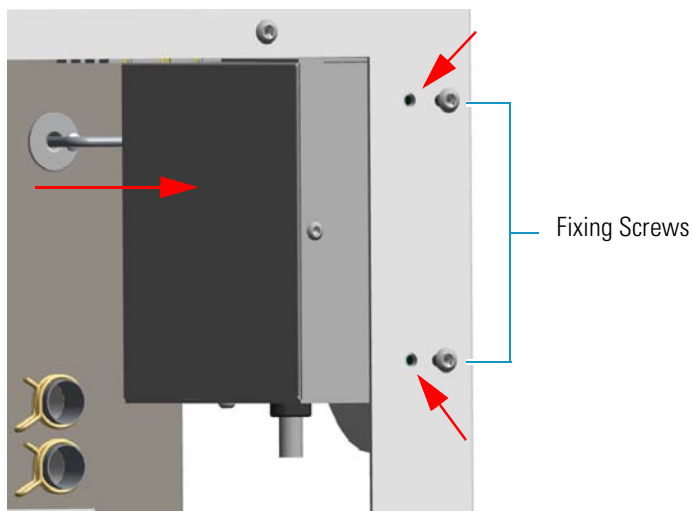


- c. Fix the hydrogen sensor using the two fixing screws. See [Figure 599](#).

**Figure 599.** Hydrogen Sensor Installation (3)

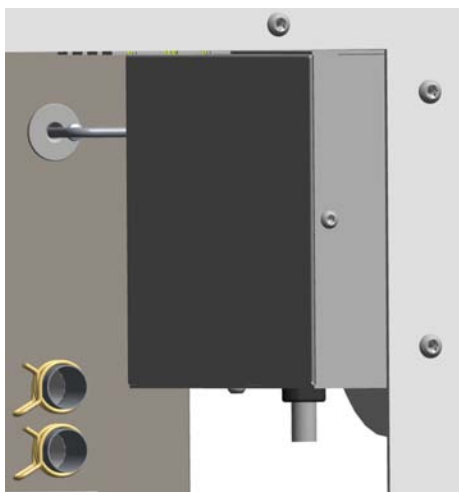

- d. Look into oven the hydrogen sensor tube that protrudes into the oven. See [Figure 600](#).

**Figure 600.** Sensor Tube Protruding Into the Oven

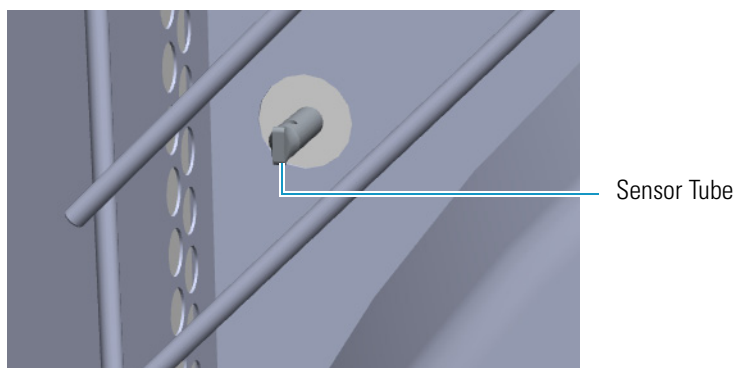



11. Connect the sensor cable to the 10-pin connector marked **J8 Hydrogen Sensor** of the Backplane board.
12. Close the oven door.
13. Reinstall the back cover.
  - a. Reconnect the ground wire to the back cover terminal.
  - b. Replace the cover proceeding in the reverse order in which it was removed.
14. Reinstall the left side panel.
  - a. Plug the ground wire to the panel.
  - b. Reinstall the panel proceeding in the reverse order in which it was removed.
15. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
16. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker), located at the back of the instrument, to the position I.
17. Set the normal injector, detector and GC working conditions.



# Upgrade Equipment

This chapter describes how to upgrade the TRACE 1300/TRACE 1310. See the *TRACE 1300 and TRACE 1310 Spare Parts Guide* for information about ordering the equipment in this chapter.

## Contents

- [Upgrading a TRACE 1300 to a TRACE 1310](#)
- [Upgrading a Stand Alone TRACE 1300/TRACE 1310 to MS Version](#)
- [Updating HMI Software From USB Stick](#)

# Upgrading a TRACE 1300 to a TRACE 1310

This section provides the instruction to upgrade a TRACE 1300 to a TRACE 1310 GC.




**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.

### ❖ To upgrade a TRACE 1300 to a TRACE 1310 GC

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.
3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
6. Remove the door cover of the TRACE 1300.
  - a. Open the front door and look for the fixing screws that secure the cover and the handle to the door. See [Figure 601](#).

**Figure 601.** Front Door Back View

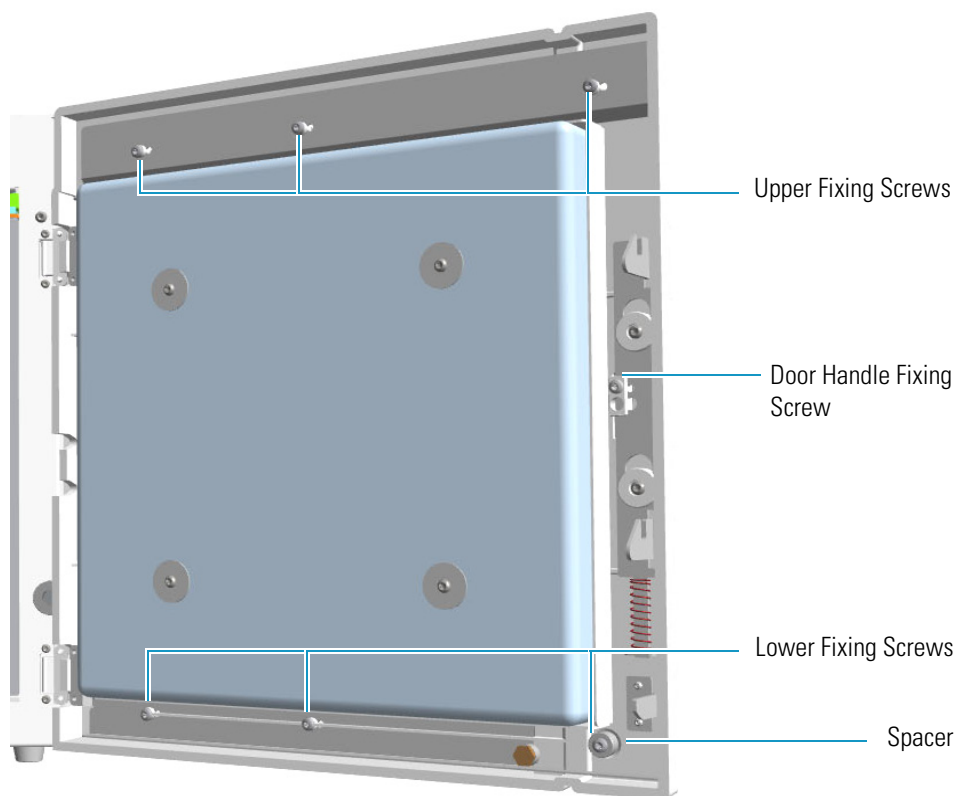


- b. Use a T20 Torxhead screwdriver to remove the screw that secure the handle to the front door.
- c. Pull the door handle out from the front door. See [Figure 602](#).  
Save the handle because it must be re-used.

## 11 Upgrade Equipment

Upgrading a TRACE 1300 to a TRACE 1310


**Figure 602.** Door Handle Removal

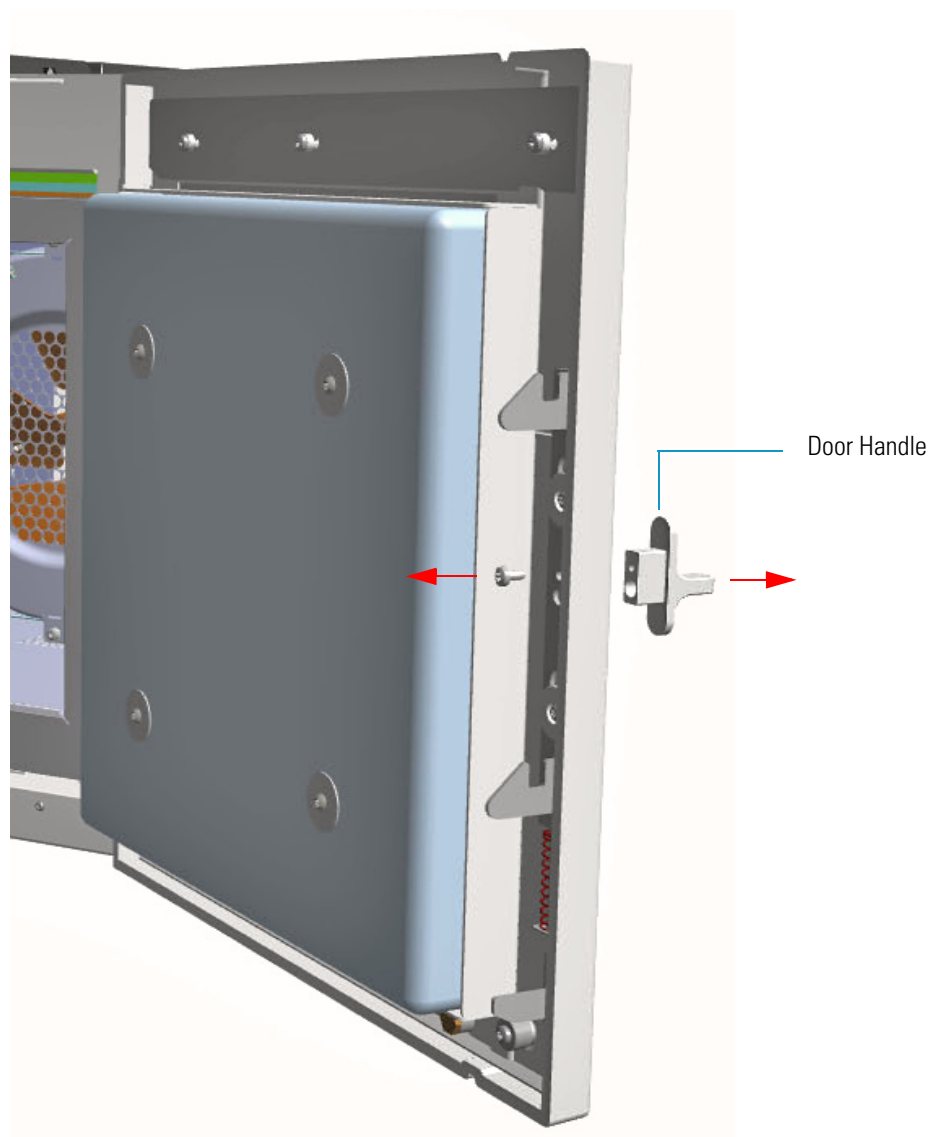

- d. Use a T20 Torxhead screwdriver to remove the three upper and the three lower screws that secure the cover to the front door. See [Figure 603](#).


**Note** The lower screw on the lower right corner is screwed into a spacer.

**Figure 603.** Door Cover Fixing Screws Removal

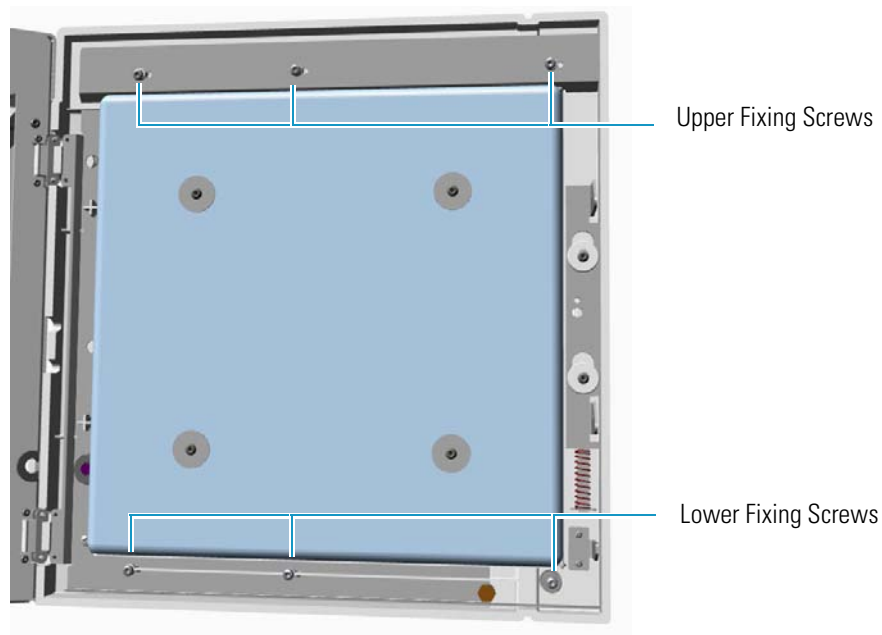

- e. Carefully pull the door cover off (see [Figure 604](#)), paying attention to the cables that connect the status panel to the internal section of the door.

**Figure 604.** Front Door Removal

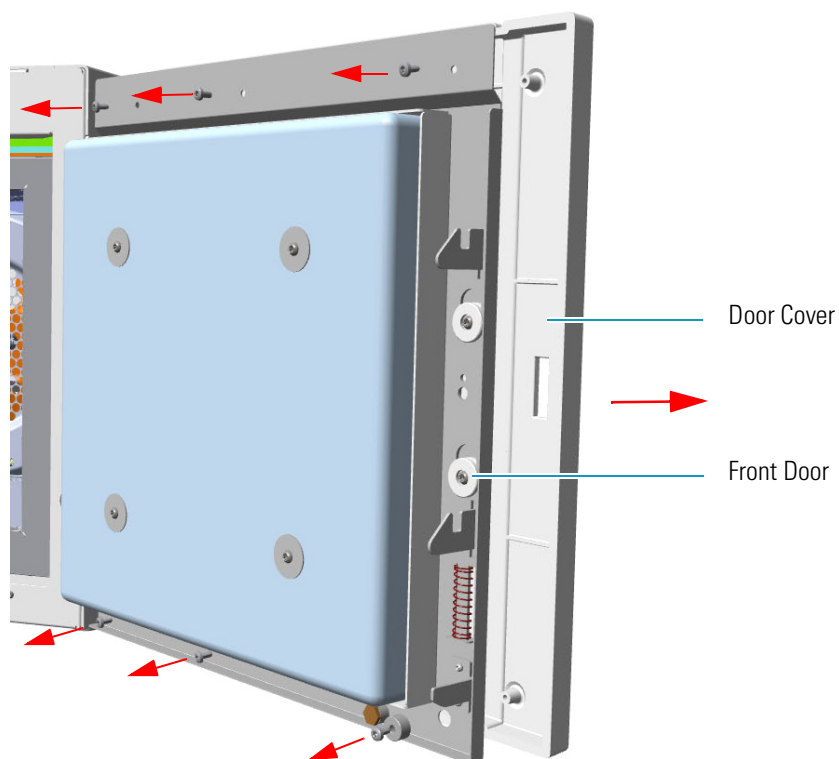


## 11 Upgrade Equipment

Upgrading a TRACE 1300 to a TRACE 1310


- f. Disconnect the flat cables from the connector located on the front of the door. See [Figure 605](#).

**Figure 605.** Cable Removal

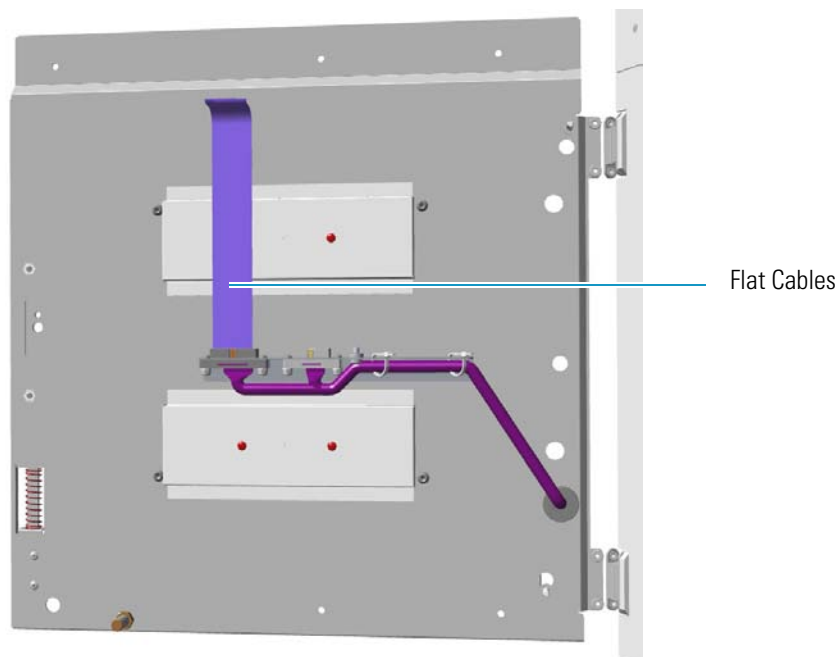

7. Install the door cover of the TRACE 1310.
  - a. The door cover of the TRACE 1310 includes the touch screen and the cables for its connection. [Figure 606](#) shows the connection points on the internal section of the door where the touch screen and ground cables must be connected.


**Figure 606.** Cables Connections Points

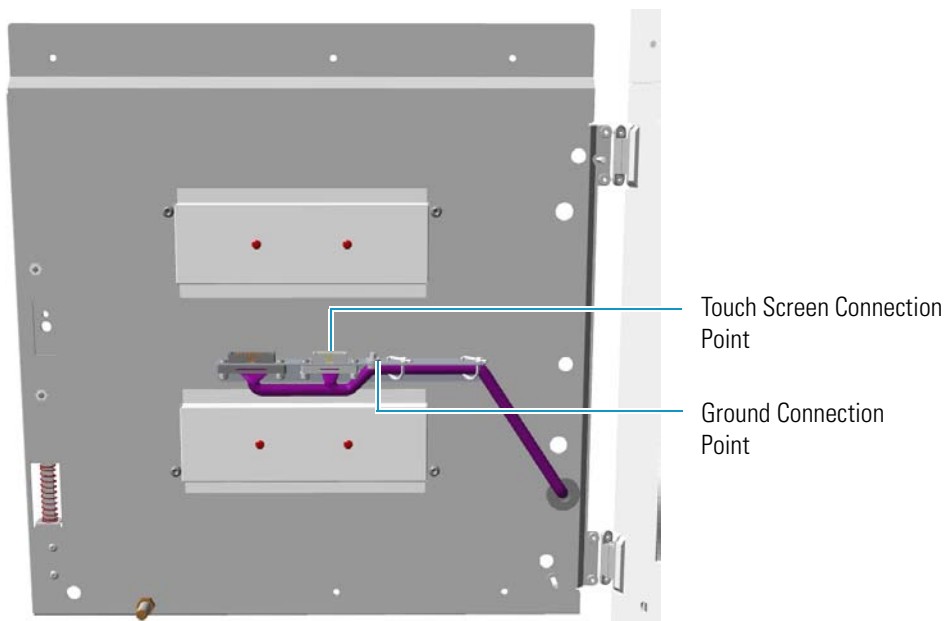



- b. Connect the touch screen and ground cables to the proper connection points located on the front of the door. See [Figure 607](#).

**Figure 607.** Cables Connections

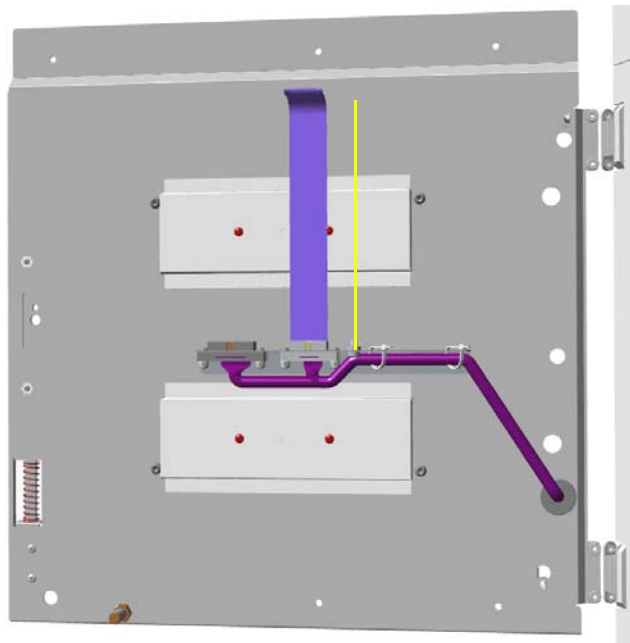

8. Mount the TRACE 1310 cover door proceeding in the reverse order in which the TRACE 1300 cover door was removed.
9. Remount the door handle.
10. Open the gas supplies.
11. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
12. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker), located at the back of the instrument, to the position I.
13. Configure the system.
14. Set the normal detector, injector, and GC working conditions.

## Upgrading a Stand Alone TRACE 1300/TRACE 1310 to MS Version

This section provides the instructions for updating your TRACE 1300/TRACE 1310 stand alone version to a MS version. According to Thermo Scientific you need a dedicated upgrade kit to couple the GC with an ISQ Series, TSQ 8000 Series, DSQ II, ITQ, or TSQ Quantum mass spectrometer.



**WARNING** This operation must be carried out by authorized and trained Thermo Fisher Scientific Service Field Engineers.

### Preliminary Operations

Before starting, the following preliminary operation must be carried out.

1. Put the GC in standby condition.
2. Cool the oven, injectors and detectors to room temperature.
3. Close the gas supplies.
4. Power off the GC.
  - a. Push down the power switch (breaker), located at the back of the instrument, to the position O.
  - b. Unplug the power cable from the AC Input connector into the back of the GC and from the wall outlet.
5. If external modules are present, unplug the power cable from the AC Input connector of each external module, and from the wall outlet.

### Getting Started

Depending on your mass spectrometer, see the following sections:

- [Coupling with the ISQ Series / TSQ 8000 Series Mass Spectrometer](#)
- [Coupling with the DSQ II Mass Spectrometer](#)
- [Coupling with the ITQ Mass Spectrometer](#)
- [Coupling with the TSQ Quantum Mass Spectrometer](#)



**IMPORTANT** The temperature control for the DSQ II, ITQ and TSQ Quantum transfer line is carried out through the **Aux Temperature/Cryo Module**. See the section [“Adding an Aux Temperature/Cryo Module”](#) on [page 443](#) for details.




**CAUTION - INSTRUMENT DAMAGE:** Condition the column before connecting it to the transfer line. The material released from the column (column bleed) during conditioning may contaminate the detector.

## Coupling with the ISQ Series / TSQ 8000 Series Mass Spectrometer

❖ **To update the GC for the coupling with the ISQ Series /TSQ 8000 Series mass spectrometer**


1. Make sure that the preliminary operations have been carried out. See [“Preliminary Operations”](#) on [page 556](#).
2. Remove the left side panel.
  - a. Open the front door of the GC.
  - b. Use a T20 Torxhead screwdriver to screw the left side panel screw from the interior front panel. See [Figure 608](#). Save the screw because it will be reused later.

**Figure 608.** Left Side Panel Fixing Screw

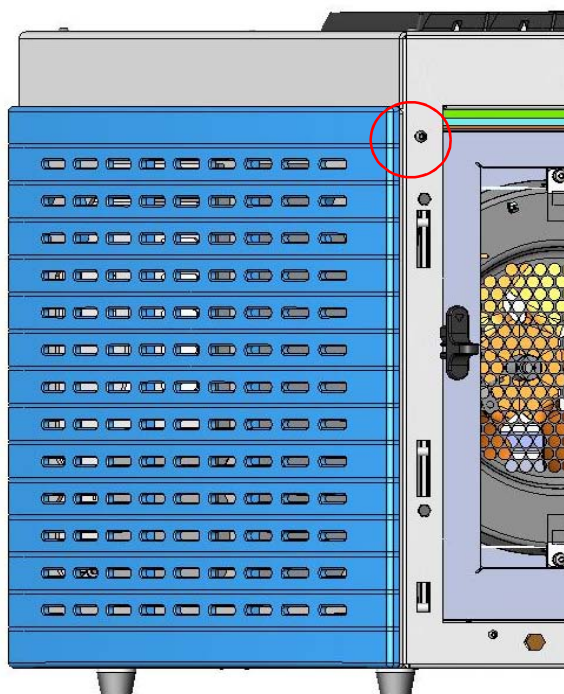


- c. Slide the panel towards the back of the instrument up to the stop.
- d. Remove the panel pulling it outward being aware that the ground wire is attached to the panel. See [Figure 609](#).
- e. Unplug the ground wire from the panel.

**Figure 609.** Left Side Panel Removal

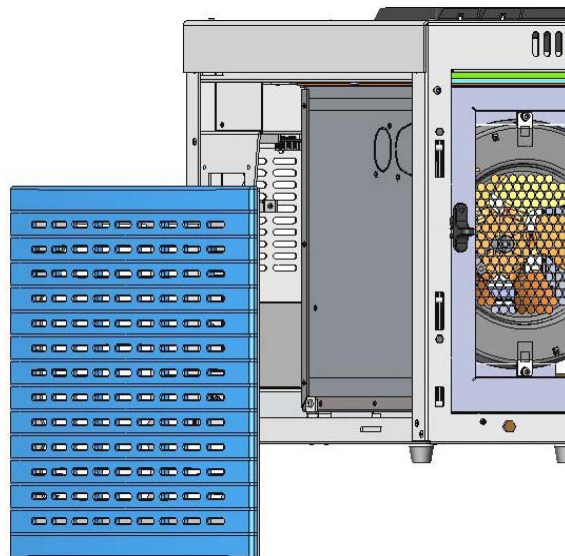


3. Prepare the duct for the transfer line inner tube.
  - a. Remove the partial cut shaped plate from the exterior wall of the oven box to access the insulating material. See [Figure 610](#).


**Figure 610.** Perform the Duct for ISQ Series / TSQ 8000 Series (1)

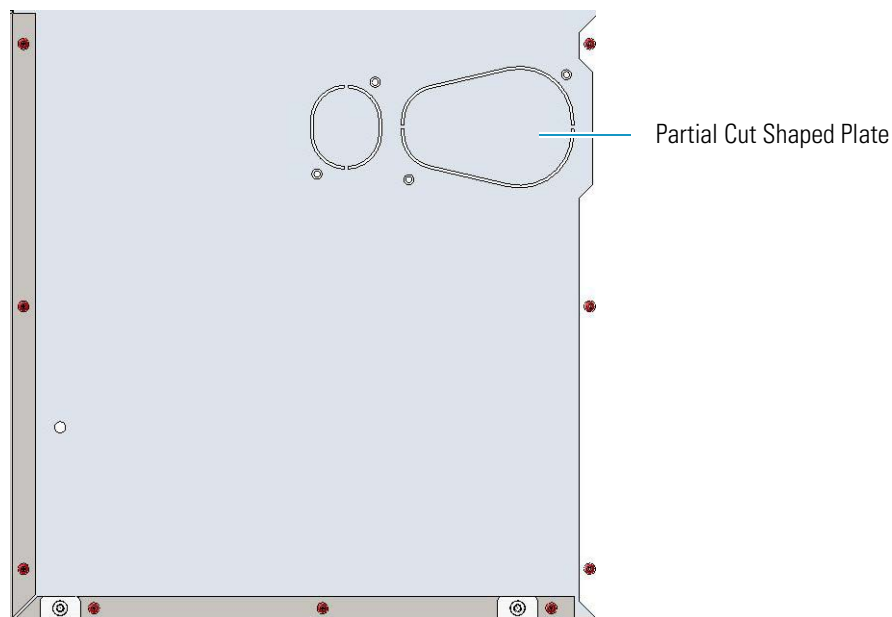

- b. Draw up the duct aligning the slot to the fixing holes on the exterior oven wall and push the duct against the insulating material up to obtain a trace. See [Figure 611](#).

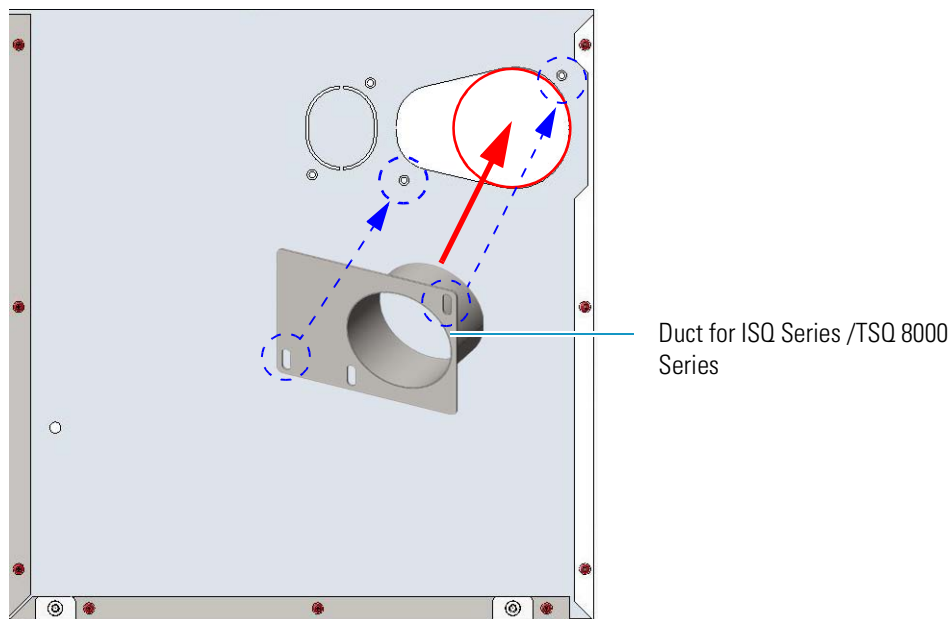
**Figure 611.** Perform the Duct for ISQ Series / TSQ 8000 Series (2)

- c. Using a knife or similar tool, gently cut the insulating material following the track.
- d. Save the removed insulating material in a safe place because it can be reused.
- e. Place the duct and fix it on the exterior wall of the oven box using the fixing screws provided. See [Figure 612](#).

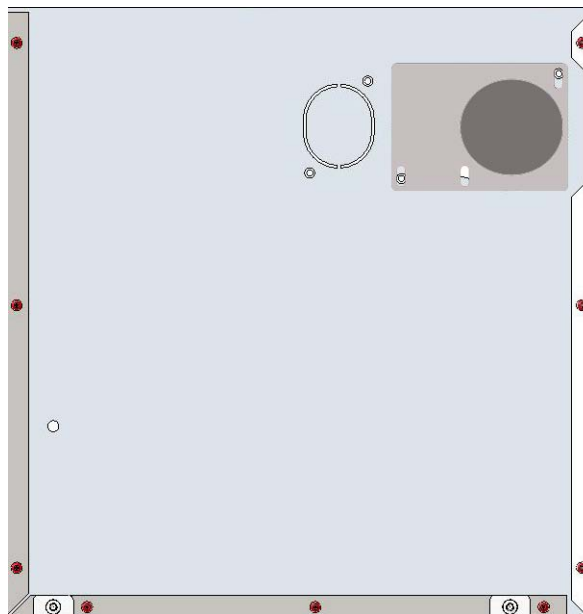**Figure 612.** Perform the Duct for ISQ Series / TSQ 8000 Series (3)


- f. On the left side wall in the interior of the oven, remove the partial cut plate from the corresponding duct. See [Figure 613](#).

**Figure 613.** Perform the Duct for ISQ Series / TSQ 8000 Series (4)

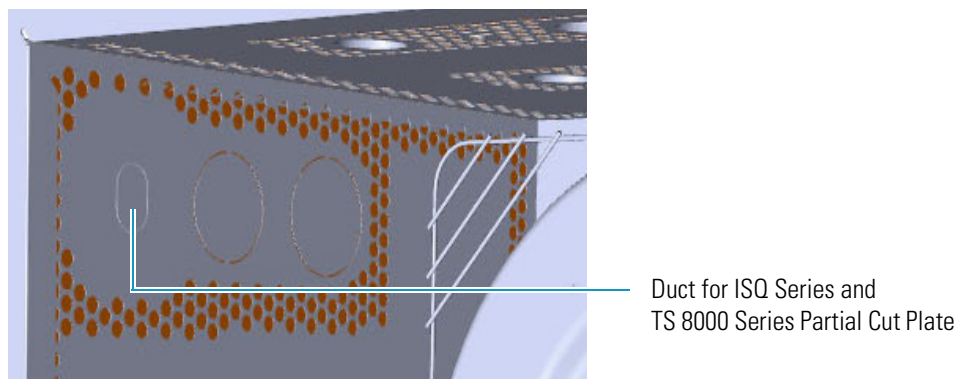

4. Replace the left side panel with the left panel for MS provided.
  - a. Remove the partial cut plate on the GC left panel for MS. See [Figure 614](#).

**Figure 614.** Left Panel for MS

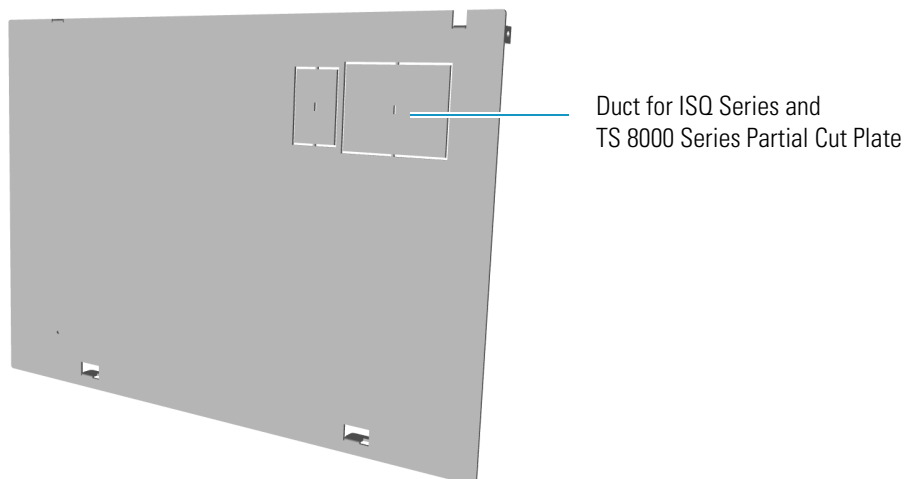

- b. Plug the ground wire previously removed to the left panel for MS provided.
  - c. Place the left panel for MS and attach the screw holding it in place.
5. Introduce the ISQ Series / TSQ 8000 Series transfer line inner tube into the oven through the duct provided.
6. Attach the transfer line to the GC column using the proper nut and ferrule.

**Note** The column must be conditioned before installing into the transfer line. See [“Installing the Column the First Time”](#) on [page 57](#)

7. Close the front door of the GC.
8. If external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
9. Power on the GC.

- a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker), located at the back of the instrument, to the position I.
10. Tune the ISQ Series / TSQ 8000 Series and set its working conditions. Set the GC working conditions accordingly.


For details please refer to the *TRACE 1300 and TRACE 1310 User Guide*, and to the *ISQ or TSQ 8000 Series User and Hardware manuals*.

## Coupling with the DSQ II Mass Spectrometer

### ❖ To update the GC for the coupling with the DSQ II mass spectrometer

1. Make sure that the preliminary operations have been carried out. See “[Preliminary Operations](#)” on [page 556](#).
2. Remove the left side panel.
  - a. Open the front door of the GC.
  - b. Use a T20 Torxhead screwdriver to screw the left side panel screw from the interior front panel. See [Figure 615](#). Save the screw because it will be reused later.

**Figure 615.** Left Side Panel Fixing Screw

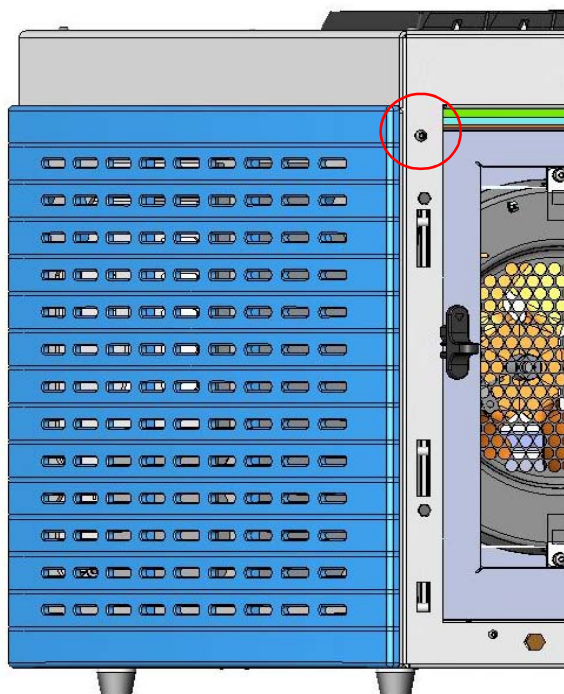


- c. Slide the panel towards the back of the instrument up to the stop.

## 11 Upgrade Equipment

Upgrading a Stand Alone TRACE 1300/TRACE 1310 to MS Version


- d. Remove the panel pulling it outward being aware that the ground wire is attached to the panel. See [Figure 616](#).
- e. Unplug the ground wire from the panel.

**Figure 616.** Left Side Panel Removal

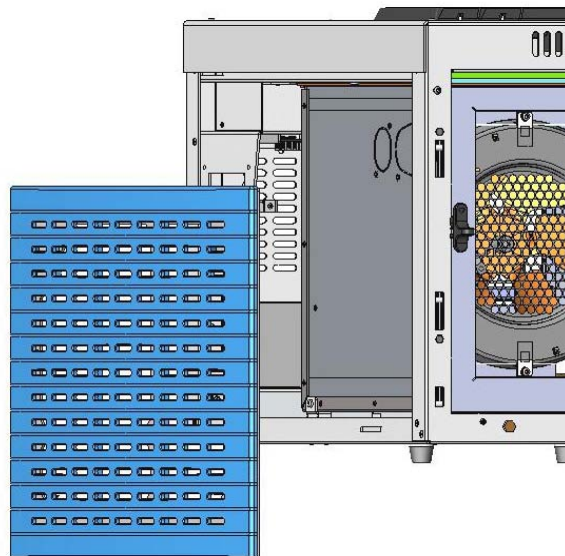

3. Prepare the duct for the transfer line inner tube.
  - a. Remove the partial cut shaped plate from the stirrer wall of the oven box to access the insulating material. See [Figure 617](#).


**Figure 617.** Perform the Duct for DSQ II (1)

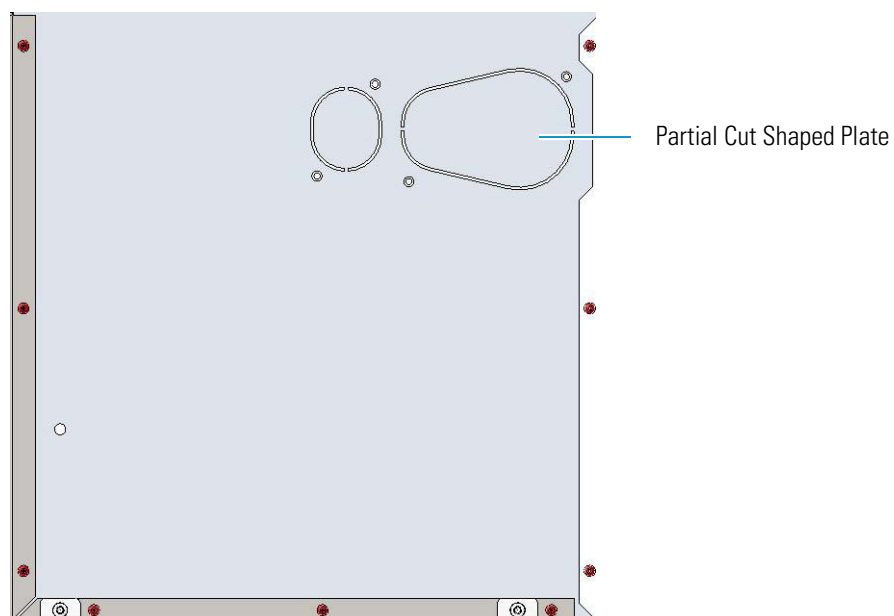




- b. Draw up the duct aligning the slot to the fixing holes on the exterior oven wall and push the duct against the insulating material up to obtain a track. See [Figure 618](#).

**Figure 618.** Perform the Duct for DSQ II (2)

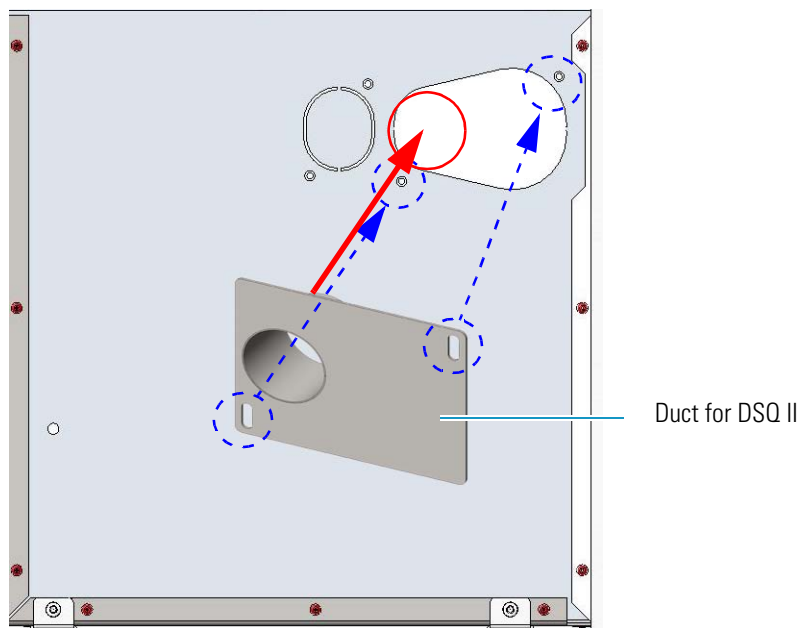


- c. Using a knife or similar tool, gently cut the insulating material following the track.
- d. Save the removed insulating material in a safe place because it can be reused.
- e. Place the duct and fix it on the exterior wall of the oven box using the fixing screws provided. See [Figure 619](#).

**Figure 619.** Perform the Duct for DSQ II (3)

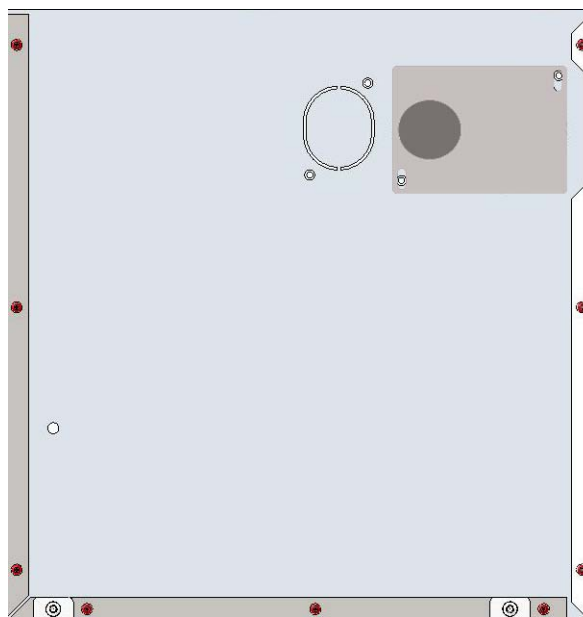

- f. On the left side wall in the interior of the oven, remove the partial cut plate from the corresponding duct. See [Figure 620](#).


**Figure 620.** Perform the Duct for DSQ II (4)

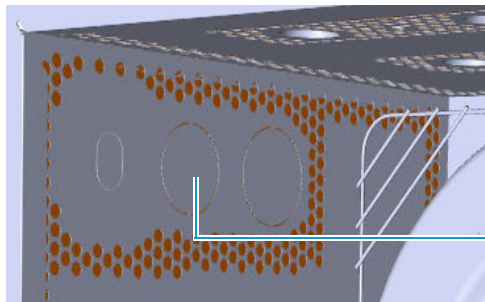

Duct for DSQ II Partial Cut Plate

4. Replace the left side panel with the left panel for MS provided.
  - a. Remove the partial cut plate on the GC left panel for MS. See [Figure 621](#).

**Figure 621.** Left Panel for MS

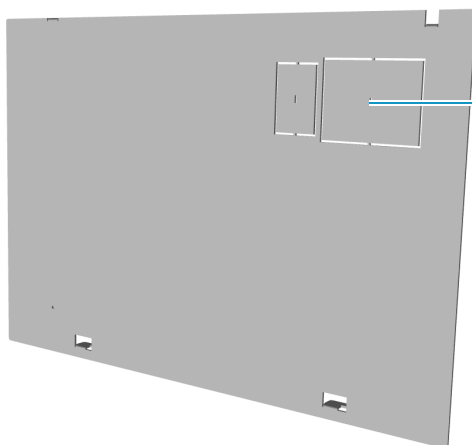

Duct for DSQ II Partial Cut Plate

- b. Plug the ground wire previously removed to the left panel for MS provided.
  - c. Place the left panel for MS and attach the screw holding it in place.
5. Introduce the DSQ II transfer line inner tube into the oven through the duct provided.
6. Attach the transfer line to the GC column using the proper nut and ferrule.


**Note** The column must be conditioned before installing into the transfer line. See [“Installing the Column the First Time”](#) on [page 57](#).

7. Place the Aux Temperature/Cryo module into the housing on the back of the GC. See the section [“Adding an Aux Temperature/Cryo Module”](#) on [page 443](#) for details.



8. Make the electrical connections. See [Figure 622](#).

**Figure 622.** Electrical Connections

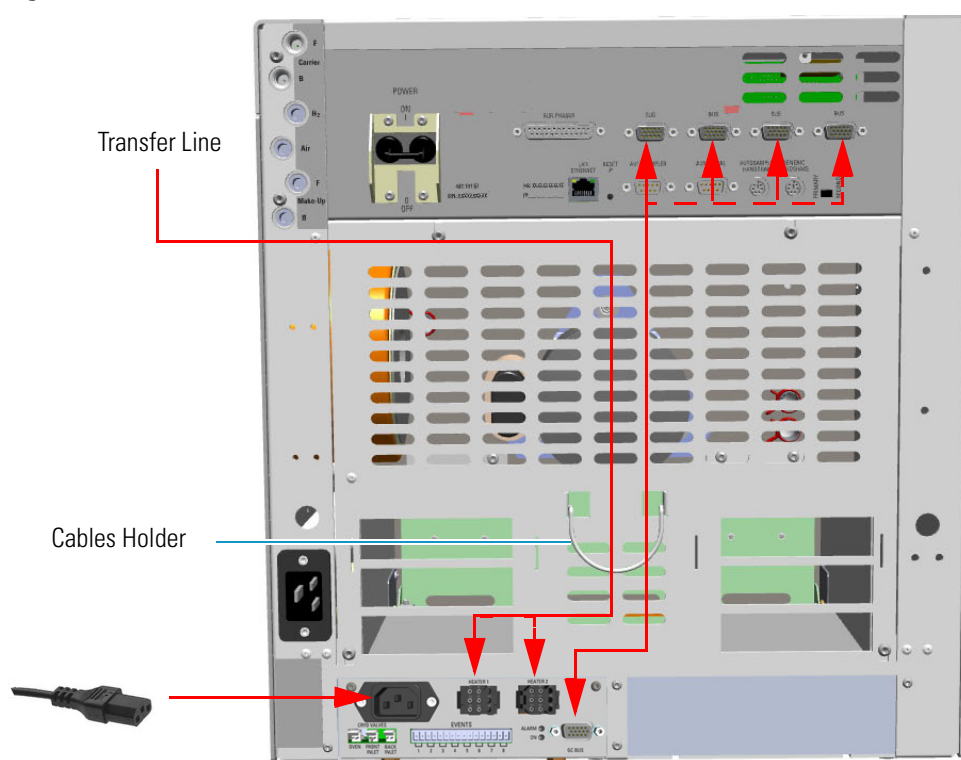

- Connect the heater cable coming from the transfer line to the connector marked **Heater 1** or **Heater 2** on the front of the module
- Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.

**Tip** To avoid the contact with the hot air from the vents, it is suggested to gather the electrical cables into the cables holder.

- Plug the power cable to the AC Input connector on the front of the module and to the wall outlet. The LED marked **On** lights.

9. Close the front door of the GC.


10. If other external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
11. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker), located at the back of the instrument, to the position I.
12. Tune the DSQ II and set its working conditions. Set the GC working conditions accordingly. For details please refer to the *TRACE 1300 and TRACE 1310 User Guide* and to the *DSQ II User and Hardware manuals*.

## Coupling with the ITQ Mass Spectrometer

### ❖ To update the GC for the coupling with the ITQ mass spectrometer

1. Make sure that the preliminary operations have been carried out. See [“Preliminary Operations”](#) on [page 556](#).
2. Remove the left side panel.
  - a. Open the front door of the GC.
  - b. Using a T20 Torxhead screwdriver, unscrew the left side panel screw from the interior front panel. See [Figure 623](#). Save the screw because it will be reused later.


**Figure 623.** Left Side Panel Fixing Screw

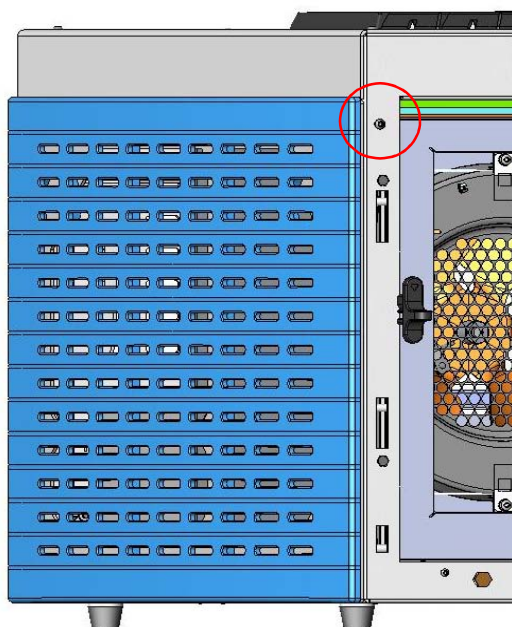

- c. Slide the panel towards the back of the instrument up to the stop.


- d. Remove the panel pulling it outward being aware that the ground wire is attached to the panel. See [Figure 624](#).
- e. Unplug the ground wire from the panel.

**Figure 624.** Left Side Panel Removal

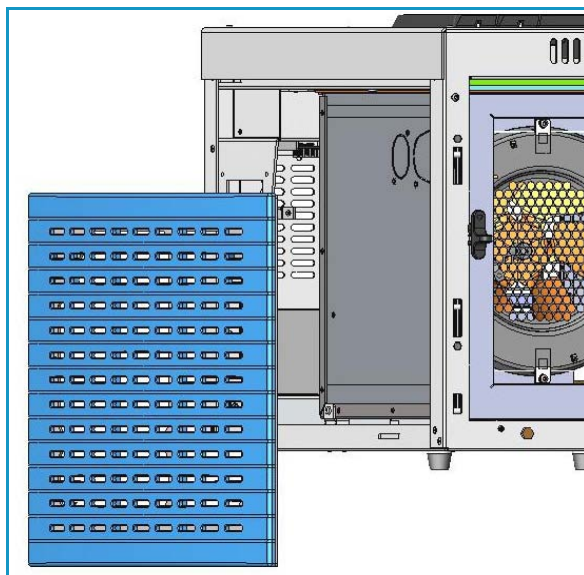


3. Prepare the duct for the transfer line inner tube.
  - a. Remove the partial cut shaped plate from the exterior wall of the oven box to access the insulating material. See [Figure 625](#).

**Figure 625.** Perform the Duct for ITQ (1)

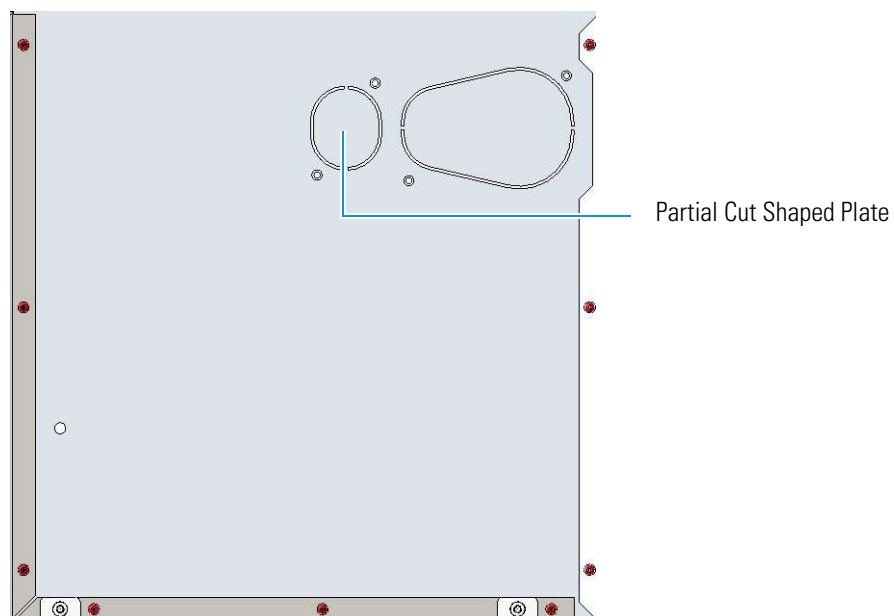


- b. Draw up the duct aligning the slot to the fixing holes on the exterior oven wall and push the duct against the insulating material up to obtain a track. See [Figure 626](#).

**Figure 626.** Perform the Duct for ITQ (2)

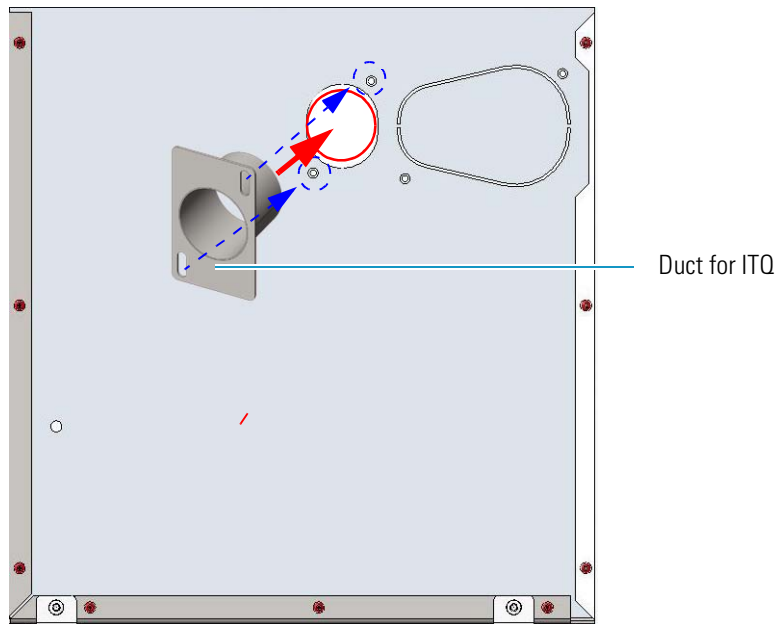


- c. Using a knife or similar tool, gently cut the insulating material following the track.
- d. Save the removed insulating material in a safe place because it can be reused.
- e. Place the duct and fix it on the exterior wall of the oven box using the fixing screws provided. See [Figure 627](#).

**Figure 627.** Perform the Duct for ITQ (3)

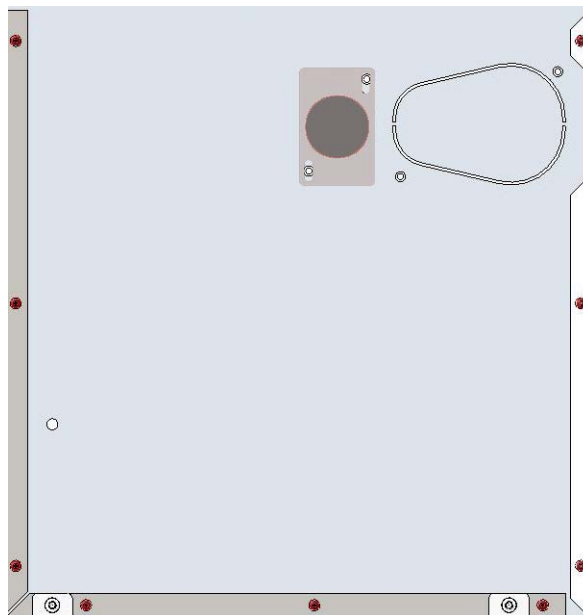


- f. On the left side wall in the interior of the oven, remove the partial cut plate from the corresponding duct. See [Figure 628](#).

**Figure 628.** Perform the Duct for ITQ (4)

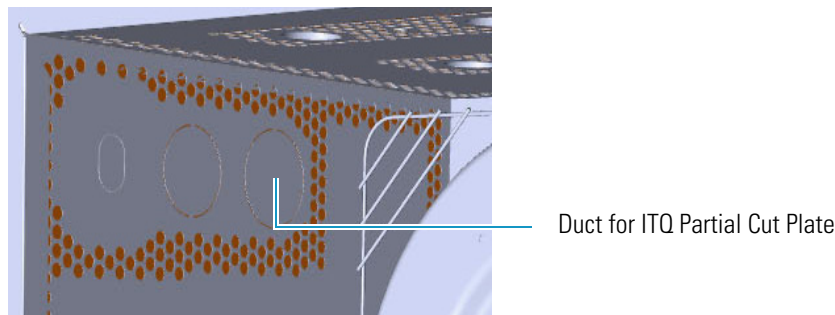

4. Replace the left side panel with the left panel for MS provided.
  - a. Remove the partial cut plate on the GC left panel for MS. See [Figure 629](#).

**Figure 629.** Left Panel for MS

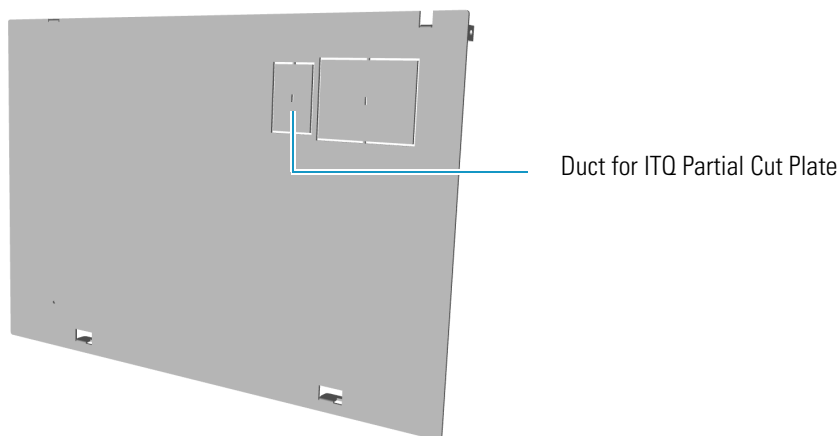


- b. Plug the ground wire previously removed to the left panel provided for MS.
  - c. Place the left panel for MS and attach the screw holding it in place.
5. Introduce the ITQ transfer line inner tube into the oven through the duct provided.
6. Attach the transfer line to the GC column using the proper nut and ferrule.

**Note** The column must be conditioned before installing into the transfer line. See [“Installing the Column the First Time”](#) on [page 57](#).

7. Place the Aux Temperature/Cryo module into the housing on the back of the GC. See the section [“Adding an Aux Temperature/Cryo Module”](#) on [page 443](#) for details.

## 11 Upgrade Equipment


Upgrading a Stand Alone TRACE 1300/TRACE 1310 to MS Version

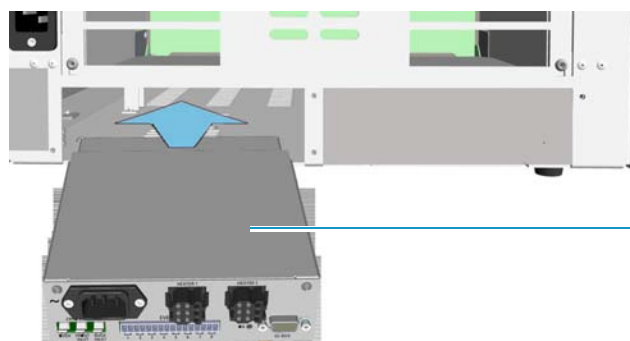

Aux Temperature/Cryo Module

8. Making the electrical connections. See [Figure 630](#).

**Figure 630.** Electrical Connections

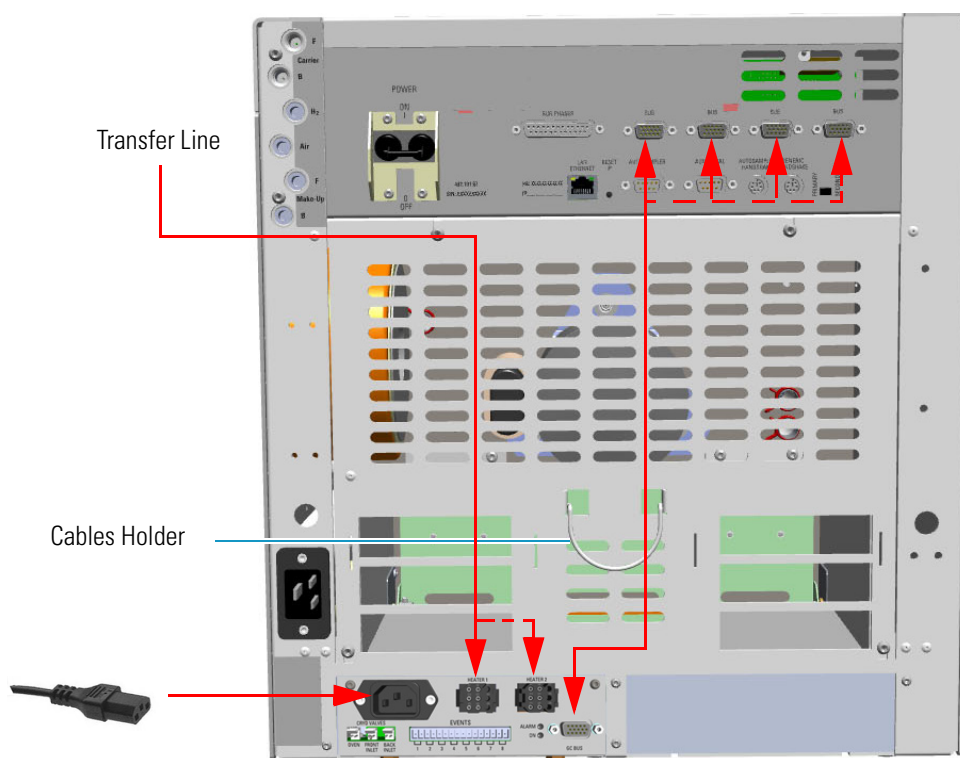

- Connect the heater cable coming from the transfer line to the connector marked **Heater 1** or **Heater 2** on the front of the module.



**ATTENTION** An extension cable transfer line is needed when connecting an ITQ mass spectrometer to the Aux Temperature/Cryo module.

- Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.

**Tip** To avoid the contact with the hot air from the vents, it is suggested to gather the electrical cables into the cables holder.




- c. Plug the power cable to the AC Input connector on the front of the module and to the wall outlet. The LED marked **On** lights.
9. Close the front door of the GC.
10. If other external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
11. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker), located at the back of the instrument, to the position I.
12. Tune the ITQ and set its working conditions. Set the GC working conditions accordingly.  
For details please refer to the *TRACE 1300 and TRACE 1310 User Guide* and to the *ITQ User and Hardware manuals*.

## Coupling with the TSQ Quantum Mass Spectrometer

### ❖ To configure the GC for the coupling with the TSQ Quantum mass spectrometer

1. Make sure that the preliminary operations have been carried out. See “[Preliminary Operations](#)” on [page 556](#).
2. Remove the right side panel.
  - a. Open the front door of the GC.
  - b. Using a T20 Torxhead screwdriver, unscrew the left side panel screw from the interior front panel. See [Figure 631](#). Save the screw because it will be reused later.

**Figure 631.** Right Side Panel Fixing Screw

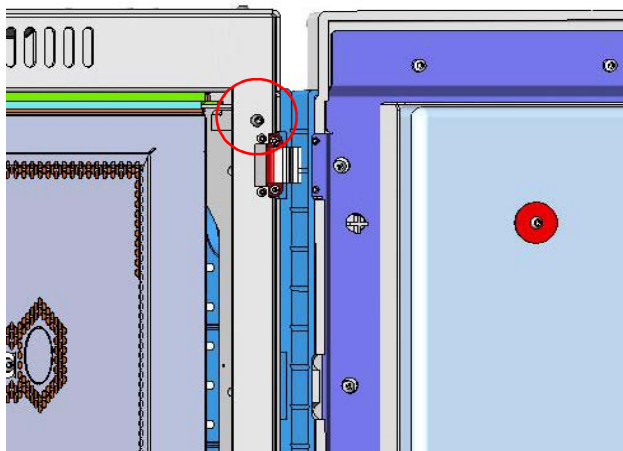


- c. Slide the panel towards the back of the instrument up to the stop.

## 11 Upgrade Equipment

Upgrading a Stand Alone TRACE 1300/TRACE 1310 to MS Version


- d. Remove the panel pulling it outward being aware that the ground wire is attached to the panel. See [Figure 632](#).
- e. Unplug the ground wire from the panel.

**Figure 632.** Right Side Panel Removal

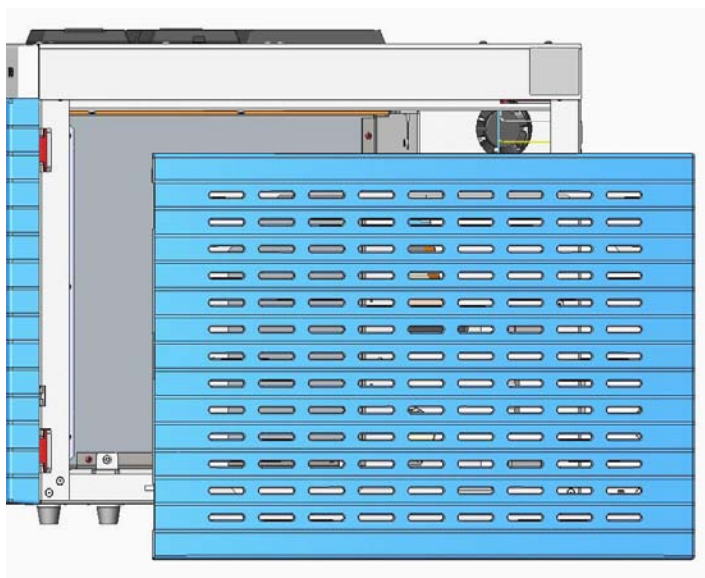


3. Prepare the duct for the transfer line inner tube.
  - a. Remove the partial cut shaped plate from the exterior wall of the oven box to access the insulating material. See [Figure 633](#).

**Figure 633.** Perform the Duct for TSQ Quantum (1)

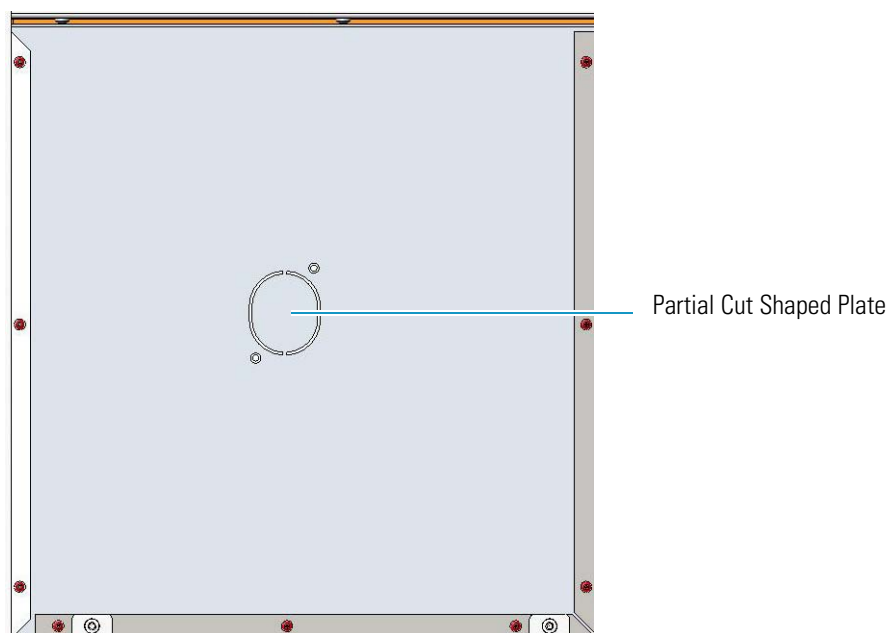


- b. Draw up the duct aligning the slot to the fixing holes on the exterior oven wall and push the duct against the insulating material up to obtain a track. See [Figure 634](#).

**Figure 634.** Perform the Duct for TSQ Quantum (2)

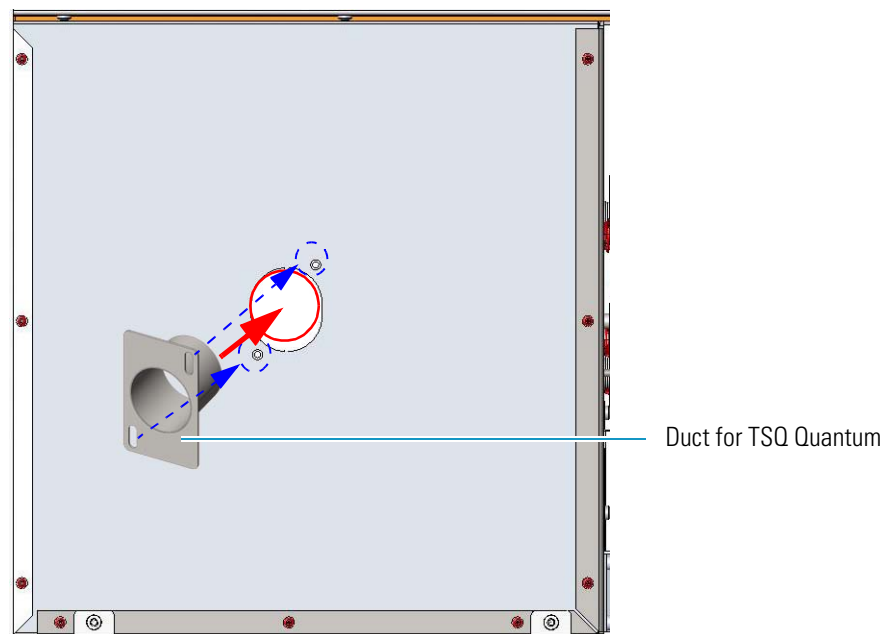


- c. Using a knife or similar tool, gently cut the insulating material following the track.
- d. Save the removed insulating material in a safe place because it can be reused.
- e. Place the duct and fix it on the exterior wall of the oven box using the fixing screws provided. See [Figure 635](#).

**Figure 635.** Perform the Duct for TSQ Quantum (3)

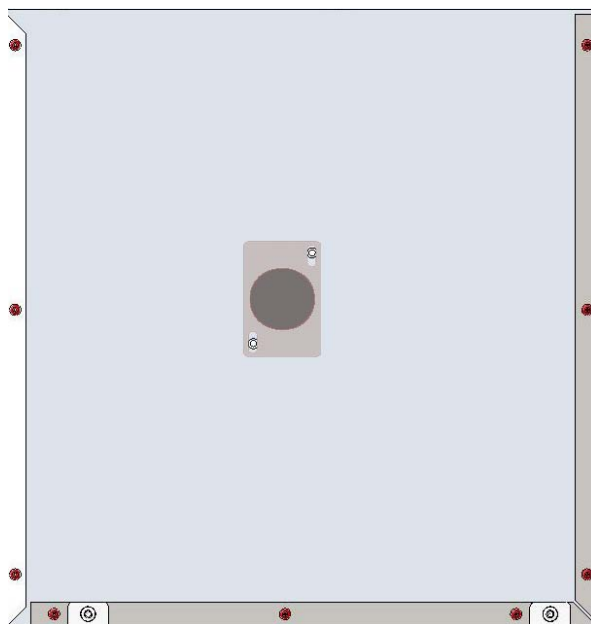


- f. On the right side wall in the interior of the oven, remove the partial cut plate from the corresponding duct. See [Figure 636](#).

**Figure 636.** Perform the Duct for TSQ Quantum (4)

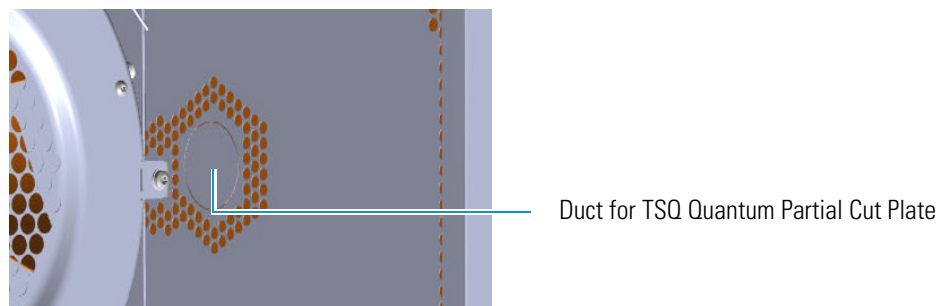

4. Replace the right side panel with the right panel provided for MS.
  - a. Using a T20 Torxhead screwdriver, unscrew the duct plate screws from the right panel for MS. See [Figure 637](#).


**Figure 637.** Right Panel for MS

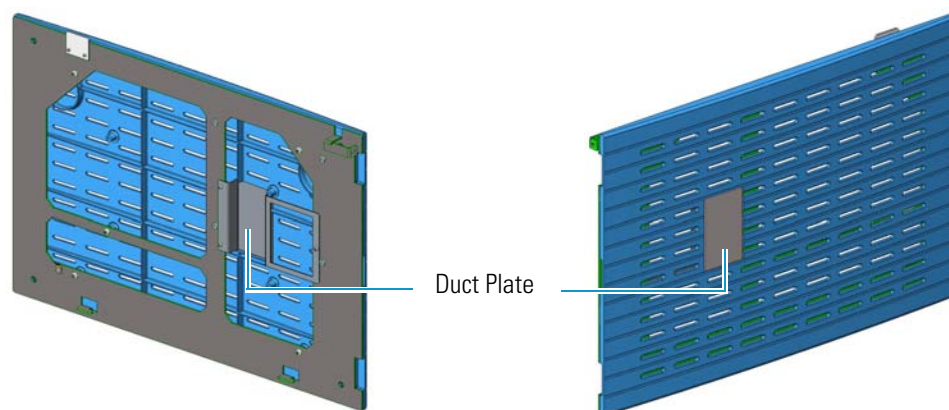

- b. Plug the ground wire previously removed to the left panel of MS provided.
  - c. Place the right panel for MS and attach the screw holding it in place.
5. Introduce the TSQ Quantum transfer line inner tube into the oven through the duct provided.
6. Attach the transfer line to the GC column using the proper nut and ferrule.


**Note** The column must be conditioned before installing into the transfer line. See [“Installing the Column the First Time”](#) on [page 57](#).

7. Place the Aux Temperature/Cryo module into the housing on the back of the GC. See the section [“Adding an Aux Temperature/Cryo Module”](#) on [page 443](#) for details.

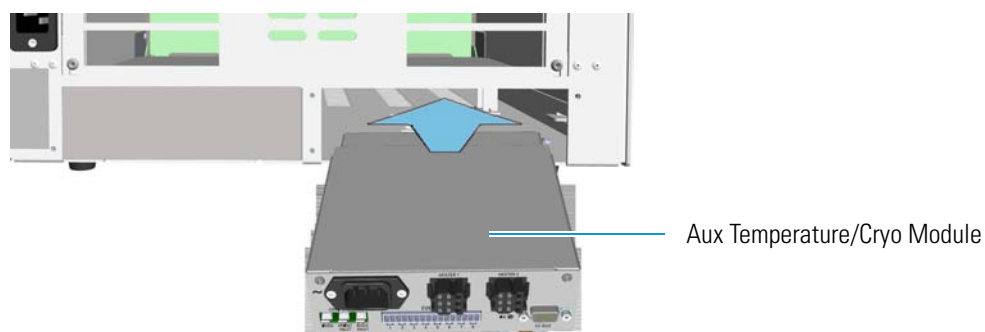

8. Making the electrical connections. See [Figure 638](#).

**Figure 638.** Electrical Connections

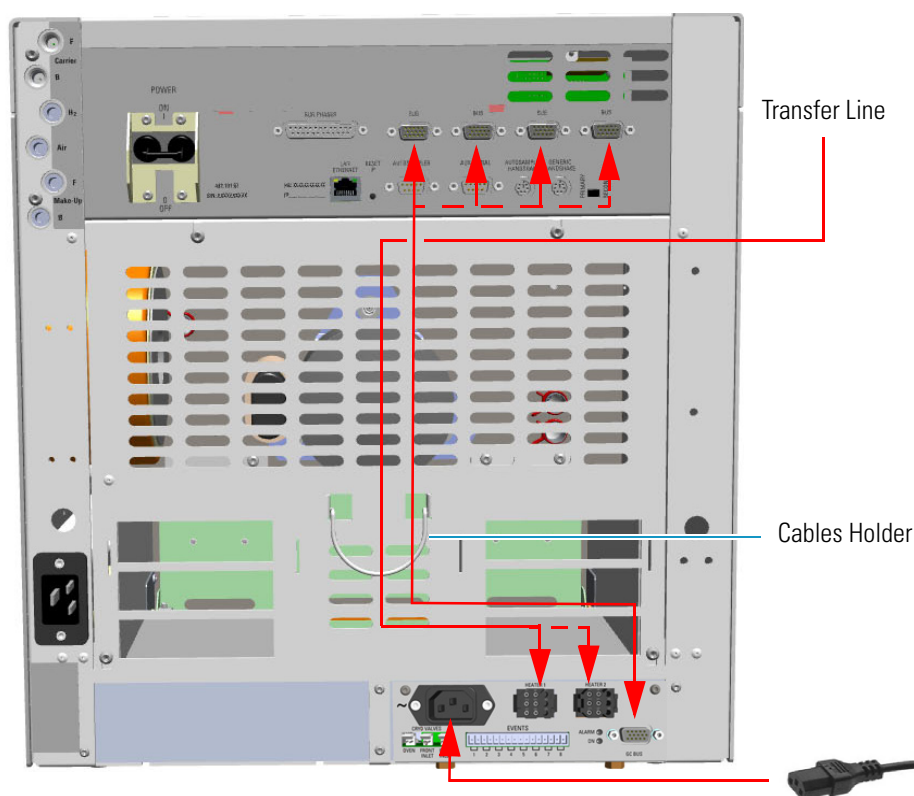

- Connect the heater cable coming from the transfer line to the connector marked **Heater 1** or **Heater 2** on the front of the module
- Using the cable provided, connect the 15-pin female connector marked **GC Bus** on the module to a **Bus** interface on the back of the GC.

**Tip** To avoid the contact with the hot air from the vents, it is suggested to gather the electrical cables into the cables holder.

- Plug the power cable to the AC Input connector on the front of the module and to the wall outlet. The LED marked **On** lights up.

9. Close the front door of the GC.
10. Open the gas supplies.
11. If other external modules are present, plug the power cable to the AC Input connector of each external module, and to the wall outlet.
12. Power on the GC.
  - a. Plug the power cable to the AC Input connector into the back of the GC and to the wall outlet.
  - b. Flip up the power switch (breaker), located at the back of the instrument, to the position I.
13. Configure and enable the cryogenic system through the user interface of your GC, or through the CDS in use. Refer to the *TRACE 1300 and TRACE 1310 User Guide*.
  - a. **Configuration and enabling through the touch screen.**
    - i. In the main menu select the **Instrument control** icon, the Instrument Control menu appears.
    - ii. In the Instrument Control menu, select the **Auxiliary** icon to open the relevant submenu.
    - iii. Set the **Aux Temp 1/2** values as required, then return to main menu.
  - b. **Configuration and enabling through the Chromatography Data System.**
    - i. In the **Configuration** window select the **Auxiliary** tab.
    - ii. Select the **Auxiliary control module** check box to enable the setting for the auxiliary control of module and option installed on your GC.
    - iii. In the **Heater 1/Heater 2** combo box specify the present heater control.
    - iv. Open the **Auxiliary setup** page.
    - v. Select **Heater 1** and/or **Heater 2** check box to enable the relevant heater control. In the adjacent field set the required temperature.
14. Set the normal detector, injector and GC working conditions.
15. Tune the TSQ Quantum and set its working conditions. Set the GC working conditions accordingly.

For details please refer to the *TRACE 1300 and TRACE 1310 User Guide* and to the *TSQ Quantum User and Hardware manuals*.

## Updating HMI Software From USB Stick

This section provides the instruction for updating the HMI software of the TRACE 1310 touch screen through an USB stick.




**ATTENTION** The upgrade must be performed by authorized and trained Thermo Fisher Scientific technical personnel.

### ❖ To update the HMI Software

1. Copy on your PC the **HMI SW XX.XX.XX.ZIP** folder of the last HMI software version received from the Thermo Fisher Scientific GC-GC/MS Custom Support.
2. Unzip the folder for extracting the files contained therein. [Figure 639](#) shows an example of the files contained into the unzipped folder.


**Figure 639.** Example of the files content into the HMI software folder

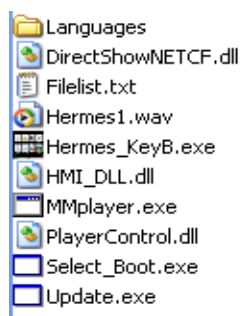

**ATTENTION** The **Hermes\_KeyB.exe** file and the **Languages** folder are the fundamental components always present in the list. The other files shown in the list might change from a software version and the next.

3. Copy **all** the files in the root of the USB stick.
4. Insert the USB stick into the USB port positioned below the touch screen.  
  
On the touch screen main menu, press the **Configuration** icon to open the Configuration menu.
5. Press the **Touch screen** icon to open the relevant menu. See [Figure 640](#).

**Figure 640.** Touch Screen Main Menu: Configuration

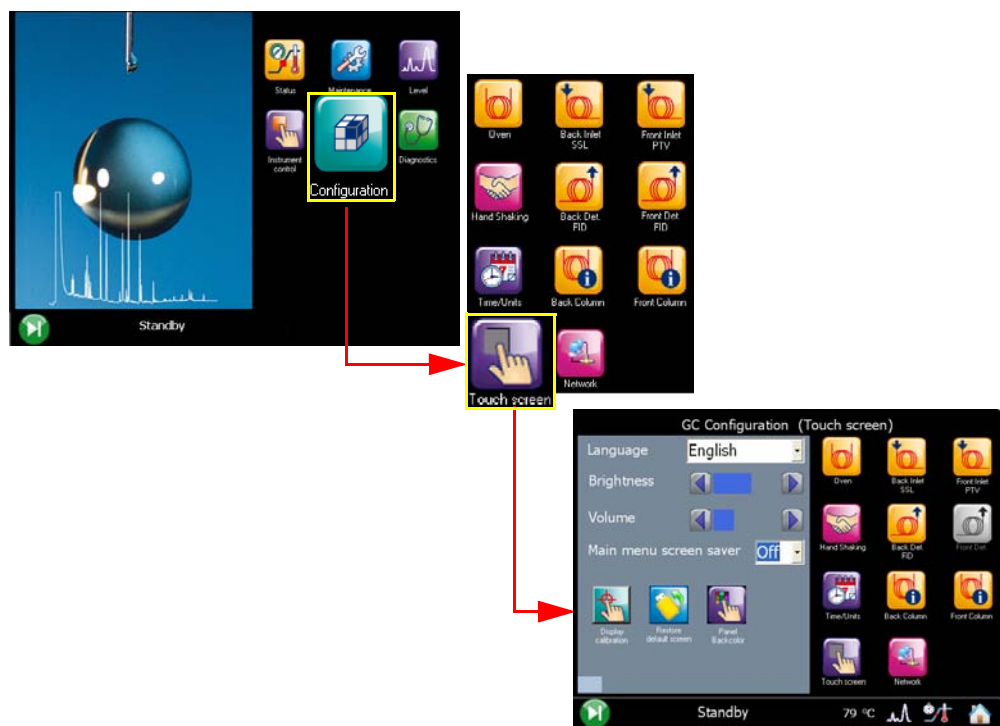

6. The icon  **Update software from USB** appears in the configuration page.  
See [Figure 641](#).

**Figure 641.** USB stick positioned into the USB port below the touch screen




**Note** It can take a few moments for USB stick to be recognized by the TRACE 1310. If you do not see the icon **Update software from USB** in the menu, then return to the Home screen, and re-enter the **Touch screen** portion of the Configuration menu.



7. Press the icon , the following page is visualized:



8. Press the icon  **Update** to start the updating process.



**WARNING** DO NOT REMOVE THE USB STICK DURING THE UPDATING PROCESS, THIS COULD DAMAGE THE INSTRUMENT.

9. At the end of the updating process, the program restarts.

10. Only now remove the USB stick from the USB port.



# Troubleshooting

In this chapter, we describe the symptoms of and remedy for each known issue with the TRACE 1300/TRACE 1310 gas chromatograph. All of these issues are related to hardware, but your instrument or software will alert you to them. For issues that you discover while reviewing your data, see the Analytical Troubleshooting section of the *TRACE 1300 and TRACE 1310 User Guide*.

## Contents

- [Investigating Power Supply Issues](#)
- [Investigating Communication Issues](#)
- [Investigating Sensitivity Issues](#)
- [Error Messages](#)
- [Contacting Technical Support](#)

## Investigating Power Supply Issues

Verify the correct power supply to the instrument.

### TRACE 1300/TRACE 1310 will not power-on

#### Possible Remedies

Make sure the TRACE 1300/TRACE 1310 power cable is properly connected to the instrument and to the correct 220/120 V main power line outlet.

Verify that the electrical outlet is functioning properly.

Power supply system is faulty. Contact your local Thermo Fisher Scientific customer support organization.

## Investigating Communication Issues

Verify the instruments is communicating with the computer.

### Software is not communicating with the TRACE 1300/TRACE 1310

#### Possible Remedies

Make sure the LAN cable is properly connected to the GC.

Confirm the TCP/IP configuration on the computer matches the GC.

Restart the GC.

### TRACE 1300/TRACE 1310 does not start or is not ready

#### Possible Remedies

Verify the GC methods and configuration.

Make sure the electrical connections have been properly carried out.

Confirm the GC handshaking parameters are set properly.

### Cannot download methods to the TRACE 1300/TRACE 1310

#### Possible Remedies

Verify that your instruments are properly configured.

## Sample data are not acquired

### Possible Remedies

Make sure the autosampler methods and configuration include starting up and injecting a sample.

You should also make sure the sample has been injected.

Make sure the cables between autosampler and GC are properly connected.

Add more disk space to the computer if necessary.

## GC is not communicating with the PC

### Possible Remedies

Make sure the GC is powered on.

Make sure the GC is properly configured.

## Autosampler is not communicating with the PC

### Possible Remedies

Make sure the autosampler is powered on.

Make sure the autosampler is properly configured.

Make sure the cable between the autosampler and PC is properly connected.

## Investigating Sensitivity Issues

Sensitivity issues are usually the result of an air leak, dirty components, or contamination. Sometimes sensitivity issues can be caused by simple problems such as the carrier gas tank running out or a sample not being injected into the GC.

If the problem is more complex, then check for air leaks or dirty components. You can prevent these problems by properly cleaning and maintaining your GC system.

It is normal to see a decrease in sensitivity in the first few injections on a clean system. Before troubleshooting for sensitivity issues, look for simple solutions, such as fixing a clogged autosampler syringe or raising the level of your sample.

## Poor sensitivity or sudden loss in sensitivity

### Possible Remedies

Check the system for leaks and address them.

Clean or replace the GC injection port liner to remove possible contamination, trim the injector end of the column, or replace the septum.

## Error Messages

Error messages are visualized in case of GC malfunctioning. See the following sections:

- [TRACE 1300 Error Messages](#)
- [TRACE 1310 Error Messages](#)

### TRACE 1300 Error Messages

**TRACE 1300** instrument malfunction, due to a component failure or to abnormal operating condition is identified by the blinking of both **Power** and **Ready** lights located on the TRACE 1300 Status Panel. See [Figure 642](#).

**Figure 642.** TRACE 1300: Alarm Identification



A typical sound (beep) is heard. Power and temperatures are automatically cut off. To reset an alarm, the GC must be powered off, and then powered on.

### TRACE 1310 Error Messages

**TRACE 1310** instrument malfunctioning, due to a component failure or to abnormal operating condition, is identified by an error messages displayed before or during the runs.



**IMPORTANT** When an alarm is displayed in the message bar, try to solve the problem, then go to the **Diagnostics** menu and press the **Reset** button to reboot the system. If the problem persists, contact the Technical Support; see “[Contacting Technical Support](#)” on [page 589](#).



- **FID Front/Back/Left/Right**
  - Unconnected
  - Temperature over limit
  - Opened PT100
  - Shorted PT100
  - Thermal safety: Not Heating or Auto Heating
  - Reset detected
- **ECD Front/Back/Left/Right**
  - Unconnected
  - PLD error
  - Temperature over limit

- Opened PT100
- Shorted PT100
- Thermal safety: Not Heating or Auto Heating
- Reset detected
- **NPD Front/Back/Left/Right**
  - Unconnected
  - Bead Current Over or Under range
  - Temperature over limit
  - Opened PT100
  - Shorted PT100
  - Thermal safety: Not Heating or Auto Heating
  - Reset detected
  - High Voltage shorted
- **TCD Front/Back/Left/Right**
  - Unconnected
  - Temperature over limit
  - Opened PT100
  - Shorted PT100
  - Thermal safety: Not Heating or Auto Heating
  - Reset detected
- **FPD Front/Back/Left/Right**
  - Unconnected
  - Temperature over limit
  - Opened PT100
  - Shorted PT100
  - Thermal safety: Not Heating or Auto Heating
  - Reset detected
  - Opened Cell PT100
  - Shorted Cell PT100
- **PDD Front/Back/Left/Right**
  - Unconnected
  - Temperature over limit

- Opened PT100
- Shorted PT100
- Reset detected
- **GDI Front/Back/Left/Right**
  - Unconnected
  - Temperature over limit
  - Opened PT100
  - Shorted PT100
  - Reset detected
- **SSL-SSLBKF Front/back**
  - Unconnected
  - Temperature over limit
  - Opened PT100
  - Shorted PT100
  - Thermal safety: Not Heating or Auto Heating
  - Reset detected
  - Loss of Carrier
- **PTV-PTVBKF Front/Back**
  - Unconnected
  - Temperature over limit
  - Opened Temperature Sensor
  - Shorted Temperature Sensor
  - Thermal safety: Not Heating or Auto Heating
  - Reset detected
  - Loss of Carrier
- **AUX Temperature EXT-V Sub Ambient**
  - Unconnected
  - X-Line B over limit
  - X-Line B Opened PT100
  - X-Line B Shorted PT100
  - X-Line A over limit
  - X-Line A Opened PT100



- X-Line A Shorted PT100
- Thermal safety: Not Heating or Auto Heating
- **Valve Oven**
  - Unconnected
  - Heater B over limit
  - Heater B Opened PT100
  - Heater B Shorted PT100
  - Heater A over limit
  - Heater A Opened PT100
  - Heater A Shorted PT100
  - Thermal safety: Not Heating or Auto Heating Heater A or B
- **AUX Carrier Primary**
  - Unconnected
  - Channel 1 Loss of Carrier
  - Channel 2 Loss of Carrier
  - Channel 3 Loss of Carrier
- **AUX Carrier Secondary**
  - Unconnected
  - Channel 1 Loss of Carrier
  - Channel 2 Loss of Carrier
  - Channel 3 Loss of Carrier
- **Oven**
  - Unconnected
  - Temperature over limit
  - Opened PT100
  - Shorted PT100
  - Thermal safety: Oven Auto Heating
  - Hydrogen Sensor Alarm
  - Hydrogen Sensor Fault
  - Reset detected
  - Thermal safety: Oven Not Heating
  - Max Temperature reached


- **Main CPU**
  - Master Safety detected
  - Main I2C Bus fault
  - Detector I2C Bus fault
  - EMI RAM Memory fault
  - Main CPU Overheating
  - Manifold Temperature Overheating
  - 48Vac missing
  - +24Vcc missing
  - +15Vcc missing
  - -15Vcc missing
  - AC Mains out of range

## Contacting Technical Support

If the information in this section does not help solve your problem, you should contact Technical Support. Be sure to reference the model, serial number, and power supply of your instrument when contacting them.


- **Serial Number** — You can find the serial number of the GC by opening the front door, and reading the serial number label on the right lower corner. See [Figure 643](#).

**Figure 643.** Identify Your Instrument: Serial Number

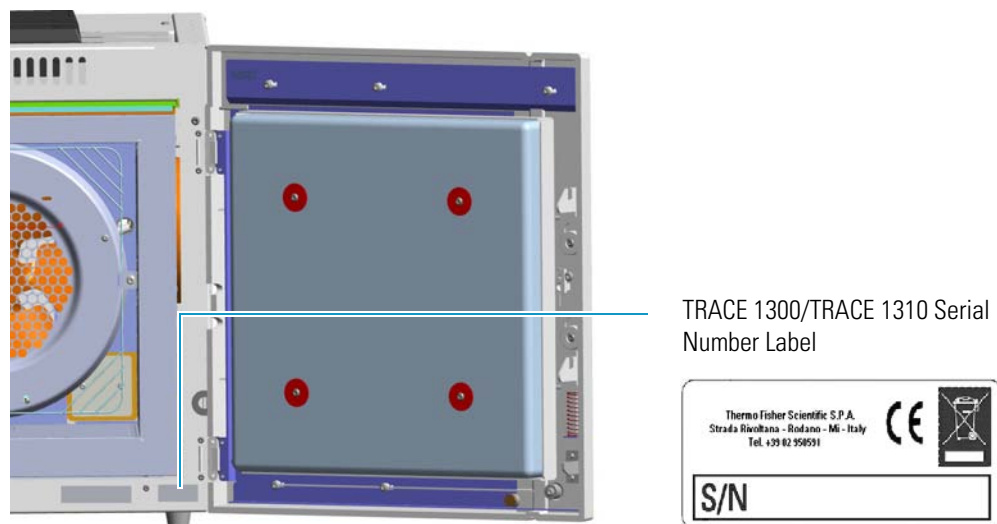

- **Power supply your GC is set** — The 120 Vac or 230 Vac power supply is indicated on the yellow label on the electronic module. See [Figure 644](#).

**Figure 644.** Identify Your Instrument: Power Supply Data

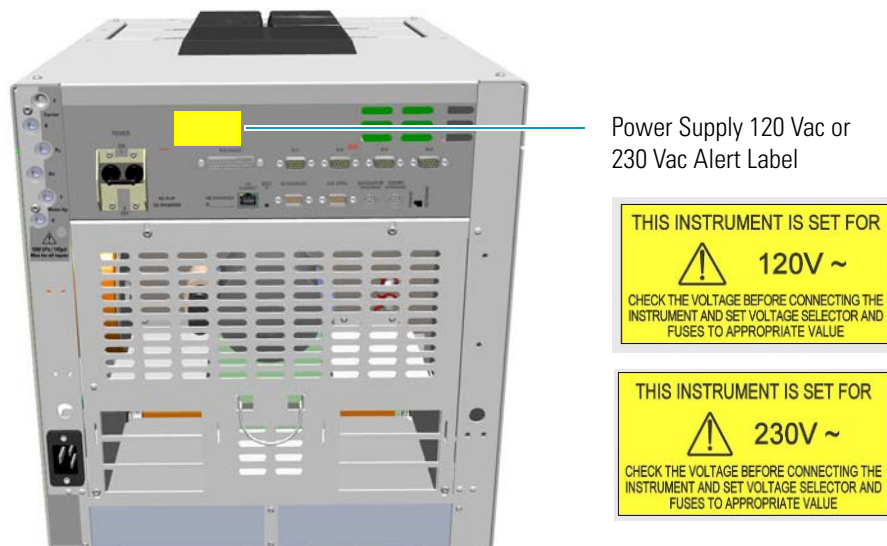

For contacting your local Thermo Fisher Scientific office or affiliate GC-GC/MS Customer Support, see the section [“Contacting Us”](#) on [page xix](#).



# Glossary

This section lists and defines terms used in this manual. It also includes acronyms, metric prefixes, and symbols.

A B C D E F G H I J K L M N O P Q R S T U V W X Y Z

## A

**A** ampere

**ac** alternating current

**ADC** analog-to-digital converter

**AOI** Analog Output Interface

## B

**b** bit

**B** byte (8 b)

**baud rate** data transmission speed in events per second

## C

**C** Carbon

**°C** Celsius

**CDS** Chromatography Data System

**CIP** Carriage and Insurance Paid To

**cm** centimeter

**CPU** central processing unit (of a computer)

**<Ctrl>** control key of the keyboard

## D

**d** depth

**DAC** digital-to-analog converter

**dc** direct current

**DS** data system

## E

**ECD** Electron Capture Detector

**EMC** electromagnetic compatibility

**ESD** electrostatic discharge

## F

**f** femto

**°F** Fahrenheit

**FID** Flame Ionization Detector

**FOB** Free on board

**FSE** Field Service Engineer

**FPD** Flame Photometric Detector

## Glossary: ft

**ft** foot

**FT-IR** fourier transform infrared spectroscopy/fourier transform infrared spectrometer

## G

**g** gram

**GC** gas chromatography- gas chromatograph

**GDI** Generic Detector Interface

**GND** electrical ground

**GSV** Gas Sampling Valve

## H

**h** height

**h** hour

**H** Hydrogen

**harmonic distortion** A high-frequency disturbance that appears as distortion of the fundamental sine wave

**He** Helium

**HeS-S/SL** Instant Connect Helium Saver Injector

**HV** high voltage

**Hz** hertz (cycles per second)

## I

**ID** inside diameter

**IEC** International Electrotechnical Commission

**Impulse** See *transient*

**in.** inch

**I/O** input/output

## K

**k** kilo ( $10^3$  or 1024)

**K** Kelvin

**kg** kilogram

**kPa** kilopascal

## L

**l** length

**L** liter

**LAN** Local Area Network

**lb** pound

**LED** light-emitting diode

## M

**m** meter (or milli [ $10^{-3}$ ])

**M** mega ( $10^6$ )

**μ** micro ( $10^{-6}$ )

**MBq** megabecquerel

**Ci** millicurie

**min** minute

**mL or ml** milliliter

**mm** millimeter

**MS** mass spectrometry-mass spectrometer

**m/z** mass-to- charge ratio

## N

**n** nano ( $10^{-9}$ )

**N** Nitrogen

**negative polarity** The inverse of a detector signal polarity.

**nm** nanometer

**NPD** Nitrogen Phosphorous Detector

## O

**OD** outside diameter

$\Omega$  ohm

## P

**p** pico ( $10^{-12}$ )

**Pa** pascal

**PCB** printed circuit board

**PDD** Pulsed Discharge Detector

**PN** part number

**psi** pounds per square inch

**PTV** Programmable Temperature Vaporizing Injector

**PTVBKF** Programmable Temperature Vaporizing  
Injector with backflush

## R

**RAM** random access memory

**<Return>** <Return> key on the keyboard

**RF** radio frequency

**ROM** read-only memory

**RS-232** industry standard for serial communication

## S

**s** second

**sag** See *surge*

**slow average** A gradual long-term change in average  
RMS voltage level, with typical duration greater than  
2 s.

**SOP** Standard Operating Procedures

**SSL** split/splitless injector

**SSLBKF** split/splitless injector with Backflush

**source current** The current needed to ignite a  
source, such as a detector lamp.

**surge** A sudden change in average RMS voltage level,  
with typical duration between 50  $\mu$ s and 2 s.

## T

**TCD** Thermal Conductivity Detector

**transient** A brief voltage surge of up to several  
thousand volts, with a duration of less than 50  $\mu$ s.

## V

**V** volt

**Vac** volts, alternating current

**Vdc** volts, direct current

**VGA** Video Graphics Array

## W

**w** width

**W** Watt

When a unit of measure has a quotient (e.g. Celsius  
degrees per minute or grams per liter) this can be written  
as negative exponent instead of the denominator:

For example:  
 $^{\circ}\text{C min}^{-1}$  instead of  $^{\circ}\text{C/min}$   
 $\text{g L}^{-1}$  instead of  $\text{g/L}$

