Bluetooth Body Scale

English

Safety precautions

Do not use the device for anything other than its intended use

Do not disassemble, modify, or repair the device


Any changes or modifications to the device can void your manufacturer's warranty. For service, take the device to a Samsung Service Center.

When using the device, mind the following:

- Humidity and liquids may damage parts or electronic circuits.
- Do not store the device in hot areas. High temperature may cause damage to the device or shorten battery life.
- Do not drop the device or cause impacts to it.
- Do not jump on the device.

- This manual's images and content may differ from the product.
- This product is only a health assisting device. It cannot be used for medical diagnostic purposes.


Overview

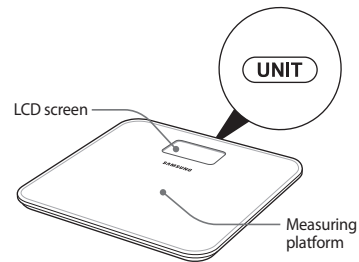Check the product box for the following items: Bluetooth Body Scale, three AAA batteries, user manual. The supplied items may vary depending on the region.

Inserting the batteries

- 1 Remove the battery cover.

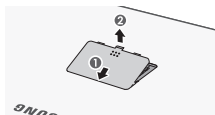- 2 Insert the batteries into the battery slots, and then replace the battery cover.

For taring weight, place the device on a flat surface within 30 seconds after inserting the batteries. If not, measured weight can be different from actual weight.

Connecting with a smartphone

- 1 Turn on the smartphone, and tap **S Health** on the Applications screen.
- 2 Tap **Health board** → **Weight** → **≡** → **Weighing accessory manager** → **EI-HS10**.
- 3 Connect with the device using the instructions on the smartphone's display.

Up to four users can be registered to the device.

Disconnecting a smartphone

- 1 Turn on the device, and press and hold **UNIT** to enter **PAIRING MODE**.
- 2 Press **UNIT** to select a user, and then press and hold **UNIT**.
- 3 Press **UNIT** to select **DISCONNECT**, and then press and hold **UNIT**.

- 4 When disconnected, the device displays **SUCCESS**.

Selecting a registered user

- 1 Step lightly on the measuring platform to turn the device on.
- 2 If the LCD screen is turned on, step on the bottom of the measuring platform to select a registered user.

When stepping on areas other than the lower part of the measuring platform without choosing a user, the device works as a standard scale. Because it acts as an unconnected scale and measurements do not transfer to the smartphone.

Measuring weight

- 1 Turn the device on, select a user, and then step on the device.

- 2 Step off the device if measuring weight is complete.
- 3 Measured weight is displayed momentarily, and is transferred to the connected smartphone.

- The measurable weight range is 10–180 kg / 22 lb–396 lb.
- To change the unit of measure (kg/lb), press **UNIT**.

FCC

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Users are not permitted to make changes or modify the device in any way.

Changes or modifications not expressly approved by Samsung will void the user's authority to operate the equipment.

Recycling Direct

To dispose electronics properly, follow the recycling process regulated in your area or go to Samsung website to find the nearest recycling location:
www.samsung.com/recyclingdirect
 Or call, (877) 278 – 0799.

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FCC STATEMENT

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

WARNING

Any changes or modifications(including the antennas) not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio / TV technician for help.

RF Exposure Statement (2.1091)

FCC RF Radiation Exposure Statement: This equipment complies with FCC RF Radiation exposure limits set forth for an uncontrolled environment. This device and its antenna must not be co-located or operating in conjunction with any other antenna or transmitter.

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IC STATEMENT

This device complies with Industry Canada licence-exempt RSS standard(s).

Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

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Samsung Telecommunications America, LLC