

TRIO™ MP-245

THREE-AXIS MOTORIZED MICROMANIPULATOR SYSTEM

WITH SYNTHETIC FOURTH "D" AXIS AND
USB INTERFACE FOR EXTERNAL CONTROL

OPERATION MANUAL

REV. 2.43 (20190130) (FW v2.4)



SUTTER INSTRUMENT®


ONE DIGITAL DRIVE
NOVATO, CA 94949

VOICE: 415-883-0128 WEB: WWW.SUTTER.COM
FAX: 415-883-0572 EMAIL: INFO@SUTTER.COM



Copyright © 2018 Sutter Instrument Company. All Rights Reserved.

TRIO™ is a trademark of Sutter Instrument Company.

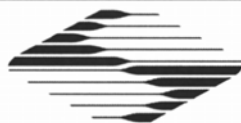

CE EU Declaration of Conformity

Application of Council Directives:
2014/30/EU (EMC), 2014/35/EU (LVD), and 2011/65/EU (RoHS 2)

Manufacturer's Name: Sutter Instrument Company

Manufacturer's Address: One Digital Drive
Novato, CA. 94949 USA
Tel: +1 415 883 0128

Equipment Tested: **TRIO MP-245** 3-Axis Motorized Micromanipulator System

Model(s): TRIO MP-245 system consisting of
TRIO-245/E (controller & ROE (Rotary Optical Encoder) for user control), with Power adapter (100-240 VAC to 24VDC) (same as, and tested as, a QUAD/E
TRIO-245/M (motorized micromanipulator electromechanical, subset of the QUAD/M),

Conforms to Standards: EMC Emissions: EN 61326-1:2013, including:
EN 55011: 2009 Group 1, Class A;
EN 61000-3-2:2015, & EN 61000-3-3:2014

EMC Immunity: EN 61000-4-2:2009, EN 61000-4-3:2011,
EN 61000-4-4:2012, EN 61000-4-5:2014,
EN 61000-4-6:2014, EN 61000-4-8:2010, &
EN 61000-4-11:2004

LVD (Safety): EN 61010-1:2010

Tested/Verified (as a QUAD system on which the TRIO MP-245 is based) by: ITC Engineering Services, Inc
9959 Calaveras Road, PO Box 543
Sunol, CA 94586-0543 USA
Tel. +1 925 862 2944 Fax: +1 925 862 9013
Email: itcemc@itcemc.com Web: www.itcemc.com
Sutter Instrument

Test Report(s): 20140120-01R1-Micromanipulator, 20140120-01
SI_EM_C_QUAD_20160713

Sutter Instrument Company hereby declares that the equipment specified above was tested and conforms to the EU Directives and Standards listed above, and further certifies conformation to the requirements of the European Union's Restriction on Hazardous Substances in Electronic Equipment Directive 2011/65/EU (RoHS 2).

Project Engineer:

Jack Belgum, Ph.D.
Senior Vice President

SUTTER INSTRUMENT®

One Digital Drive, Novato, CA 94949 USA Phone: +1 415 883 0128 Fax: +1 415 883 0572
Email: info@sutter.com Web: <http://www.sutter.com>

DISCLAIMER

The **TRIO MP-245** consists of one electromechanical micromanipulator device and one ROE (Rotary Optical Encoder) with integrated controller. The purpose of the system is for the manipulation at the micro level of micropipettes and probes used in conjunction with a microscope. No other use is recommended.

This instrument is designed for use in a laboratory environment. It is not intended nor should it be used in human experimentation or applied to humans in any way. This is not a medical device.






Unless otherwise indicated in this manual or by Sutter Instrument Technical Support for reconfiguration, do not open or attempt to repair the instrument.

Do not allow an unauthorized and/or untrained operative to use this device.


Any misuse will be the sole responsibility of the user/owner and Sutter Instrument Company assumes no implied or inferred liability for direct or consequential damages from this instrument if it is operated or used in any way other than for which it is designed.

SAFETY WARNINGS AND PRECAUTIONS

Electrical

- Operate the TRIO MP-245 using 110 – 240 VAC., 50-60 Hz line voltage. This instrument is designed for use in a laboratory environment that has low electrical noise and mechanical vibration. Surge suppression is recommended at all times
-  NOTE: There are no user-replaceable fuses in the TRIO MP-245 system.
-   The TRIO MP-245 system's power supply consists of an external AC to DC switching power adapter. If the external power adapter is damaged due to a mains over or under voltage, it must be replaced.
-   **GROUNDING/EARTHING:** Proper grounding protects the ROE/controller electronics, reduces/eliminates electromagnetic interference, and improves the safety of the system operator. The ROE/controller provides a socket (labeled GROUND) that accepts a banana plug attached to a suitably gauged insulated wire, the other end of which (alligator clip) connects to a solid, proper ground.

Avoiding Electrical Shock and Fire-related Injury

- Always use the grounded power cord provided to connect the system's power adapter to a grounded/earthed mains outlet (3-prong). This is required to protect you from injury in the event that an electrical hazard occurs.
- Do not disassemble the system. Refer servicing to qualified personnel.
-  To prevent fire or shock hazard do not expose the unit to rain or moisture.

Electromagnetic Interference

To comply with FDA and CE/EU electromagnetic immunity and interference standards; and to reduce the electromagnetic coupling between this and other equipment in your lab always use the type and length of interconnect cables provided for interconnecting the electro-mechanical devices and ROE/controller (refer to Technical Specifications for more details).

Operational

Failure to comply with any of the following precautions may damage this device.

- This instrument is designed for operation in a laboratory environment (Pollution Degree I) that is free from mechanical vibrations, electrical noise and transients.
- **⚠ DO NOT CONNECT OR DISCONNECT THE CABLES BETWEEN THE CONTROLLER AND THE MECHANICAL UNITS WHILE POWER IS ON.**
Please allow at least 20 seconds after turning the unit off before disconnecting the mechanical units. Failure to do so may result in damage to the electronics.
- Operate this instrument only according to the instructions included in this manual.
- Do not operate if there is any obvious damage to any part of the instrument.
- **⚠ Do not operate this instrument near flammable materials.** The use of any hazardous materials with this instrument is not recommended and, if undertaken, is done so at the users' own risk.
- **⚠ Do not operate if there is any obvious damage to any part of the instrument.** Do not attempt to operate the instrument with the TRIO MP-245/M electromechanical manipulator shipping tape in place or severe motor damage may result. When transporting the mechanical manipulator, be sure to reinstall the shipping tape (using masking tape or equivalent only) to the original locations. Failure to do this may result in damage to the motors.
- **⚠ Never touch any part of the micromanipulator electromechanical device while it is in operation and moving.** Doing so can result in physical injury (e.g., fingers can be caught and pinched between the moving parts of the micromanipulator).
- **⚠ If the TRIO MP-245 system is used in a microinjection environment, please observe the following.** As with most micromanipulation devices, sharp micropipettes can fly out of their holder unexpectedly. Always take precautions to prevent this from happening. Never loosen the micropipette holder chuck when the tubing is pressurized, and never point micropipette holders at yourself or others. Always wear safety glasses when using sharp glass micropipettes with pressure tubing.
- **⚠ Take care to ensure no cables pass close to the TRIO MP-245/M electromechanical micromanipulator within the spherical movement limits of all its axes combined.**

Other

- Retain the original packaging for future transport of the instrument.
- Sutter Instrument reserves the right to change specifications without prior notice.
- Use of this instrument is for research purposes only.

Handling Micropipettes

⚠ Failure to comply with any of the following precautions may result in injury to the users of this device as well as those working in the general area near the device.

- The micropipettes used with this instrument are very sharp and relatively fragile. Avoid contact with micropipette tips to prevent accidentally impaling oneself.
- Always dispose of micropipettes by placing them into a well-marked, spill-proof "sharps" container.

TABLE OF CONTENTS

DISCLAIMER	3
SAFETY WARNINGS AND PRECAUTIONS	3
Electrical	3
Avoiding Electrical Shock and Fire-related Injury	3
Electromagnetic Interference	3
Operational	3
Other	4
Handling Micropipettes.....	4
1. INTRODUCTION	9
1.1 Structure of the TRIO MP-245 Documentation Package	9
1.2 Components of the TRIO MP-245 System.....	9
1.3 Overview	10
1.3.1 Features.....	10
1.3.2 Description	10
2. INSTALLATION	11
2.1 Mounting Instructions	11
2.1.1 Mounting the TRIO MP-245/M to the Stand or Platform	11
2.2 Headstage Mounting	13
2.3 Other Accessories.....	13
2.4 Electrical Connections and Initial Operating Instructions	13
2.5 ROE/Controller Rear Panel Controls and Configuration.....	14
2.5.1 Power Switch	14
2.5.2 Configuration Switches	15
2.5.2.1 Switches 1, 2, 3 and 4	15
2.5.2.2 Switch 5: Y-Axis Lockout during Homing	15
2.5.2.3 Switch 6: Sensor Test (Firmware < v2.2)	15
2.5.2.4 Switch 6: Sensor Test (Firmware < v2.4)	16
2.5.2.5 Switch 6: Calibration Homing on Power On (Firmware v2.4+).....	16
2.5.2.6 Switches 7 through 10.....	16
2.5.2.7 Switch 7: PULSE Button Functionality (Firmware v2.4+)	16
2.5.2.8 Switch 8: Electromechanical Device Whole Step Size (Firmware v2.4+)	17
2.5.2.9 Switch 9: Y-Axis Travel Length (Firmware v2.4+).....	17
2.5.2.10 Switch 10: Length of X Axis (Firmware v2.4+)	17
3. OPERATIONS	18
3.1 Main Controls and Indicators on the ROE/Controller	18
3.2 Display.....	19
3.2.1 Initial Startup	19
3.3 Control Operations	19
3.3.1 Maximum Positive Position Values:	19
3.3.2 Setting Position for HOME or WORK.....	20
3.3.3 Setting the Angle of the Pipette/Headstage Holder	20
3.3.4 Operating the Virtual D Axis	20
3.3.5 Moving to the Home Position	20

3.3.6 Moving to the Work Position	21
3.3.7 Setting Absolute/Relative Coordinates Mode.....	21
3.3.8 Mode Indications	22
3.3.9 Speed Control and ROE Knob Movements (SPEED).....	22
3.3.10 Movement Knobs Disabling and Lock Mode ([SPEED]/LOCK)	22
3.3.11 Pausing Home Movements (HOME (while moving to Home))	22
3.3.12 Pausing Work Movements (WORK (while moving to Work))	22
3.3.13 Pulse Mode and Virtual D-Axis Movement (PULSE).....	23
3.4 Micropipette/Headstage Exchange	23
4. EXTERNAL CONTROL.....	23
4.1 General.....	23
4.2 Virtual COM Port (VCP) Serial Port Settings.....	24
4.3 Protocol and Handshaking	24
4.4 Command Sequence Formatting	24
4.5 Axis Position Command Parameters.....	25
4.6 Microsteps and Microns (Micrometers).....	25
4.7 Commands	25
4.7.1 Get Current Position ('c' or 'C') Command	25
4.7.2 Move to Controller-Defined HOME Position ('h') Command.....	26
4.7.3 Move to Controller-Defined WORK Position ('w') Command	26
4.7.4 Move to Specified "Home" Position ('H') Command	27
4.7.5 Move to Specified "Work" Position ('W') Command.....	27
4.7.6 Move in Straight Line to Specified Position at Specified Speed ('S') Command	28
4.7.7 Interrupt Straight-Line Move ('^C') Command	28
4.7.8 Move to Specified X-Axis Position ('x' or 'X') Command.....	28
4.7.9 Move to Specified Y-Axis Position ('y' or 'Y') Command.....	29
4.7.10 Move to Specified Z-Axis Position ('z' or 'Z') Command	29
4.7.11 Setting the Angle ('A') Command	29
4.7.12 Notes.....	29
5. MAINTENANCE.....	31
6. RECONFIGURATION	31
6.1 Changing the Rotary Knob Functions on the ROE/Controller	31
APPENDIX A. LIMITED WARRANTY	32
APPENDIX B. ACCESSORIES.....	33
APPENDIX C. TECHNICAL SPECIFICATIONS	34
APPENDIX D. QUICK REFERENCE	36
D.1. Manual Operation.....	36
D.2. Configuration.....	36
D.3. External Control	37
INDEX	41

TABLE OF FIGURES

Figure 1-1. The TRIO MP-245 system.....	9
Figure 2-1. Side view of TRIO MP-245/M showing mounting adapter plate and lock screws.....	12
Figure 2-2. Mounting the TRIO MP-245/M on the Adapter Plate.....	12
Figure 2-3. Rear of TRIO MP-245 ROE/Controller cabinet	14
Figure 2-4. Configuration switches on rear of TRIO MP-245 ROE/Controller unit (switch positions shown are factory defaults).....	15
Figure 3-1. LCD Display showing startup screen.	19
Figure 3-2. Startup screen	19
Figure 3-3. Factory default startup (Home) position.....	19
Figure 3-4. Maximum positive values.....	19
Figure 3-5. Moving to Home position (screen is amber while moving)	20
Figure 3-6. Factory default Home position.....	20
Figure 3-7. Example Home position defined and saved	20
Figure 3-8. Example Work position	21
Figure 3-9. Relative mode.....	21
Figure 3-10. Relative mode.....	21
Figure 3-11. Absolute mode.....	22
Figure 3-12. Angled side view of TRIO MP-245/M to change headstage mount.....	23
Figure 6-1. Locations of the axis connectors inside the ROE/Controller	31

TABLE OF TABLES

Table 2-1. Configuration Switches 1 - 4: Configuring the direction of each axis.....	15
Table 2-2. Configuring the Homing Y-Movement Lock Out.	15
Table 2-3. Configuring the Sensor Test (Firmware < v2.2).	16
Table 2-4. Configuration Switch 6: Configuring the Sensor Test.	16
Table 2-5. Configuration Switch 6: Calibration Homing on Power On (Firmware v2.4+).	16
Table 2-6. Configuration Switch 7: PULSE button functionality.....	17
Table 2-7. Configuration Switch 8: Electromechanical device whole step size (1.5 vs. 2 μ m).	17
Table 2-8. Configuration Switch 9: Y-Axis travel length (Firmware v2.4+).	17
Table 2-9. Configuration Switch 10: Length of X Axis (Firmware v2.4+).....	18
Table 3-1. Maximum positive position value of each axis	19

Table 3-2. Screen colors and modes	22
Table 4-1. USB-VCP interface serial port settings.	24
Table 4-2. Microns/microsteps conversion.	25
Table 4-3. Ranges (travel distances and bounds of each axis for compatible devices).	25
Table 4-4. Get Current Position ('c' or 'C') command.	26
Table 4-5. Move to controller-defined HOME position ('h') command.	26
Table 4-6. Move to controller -defined WORK position ('w') command.	26
Table 4-7. Move to specified "Home" position ('H') command.	27
Table 4-8. Move to specified "Work" position ('W') command.	27
Table 4-9. Straight-line move to specified position ('S') command.	28
Table 4-10. Interrupt a straight-line move in progress (' ^ C') command.	28
Table 4-11. Move to specified X-axis position ('x' or 'X') command.	28
Table 4-12. Move to specified Y-axis position ('y' or 'Y') command.	29
Table 4-13. Move to specified Z-axis position ('z') command.	29
Table 4-14. Set the angle ('A') command.	29
Table C-1. TRIO MP-245 cables and receptacles/connectors.	35
Table D-1. Configuration Switches 1 – 5.	36
Table D-2. Configuration Switches 6 – 10 (Ver. <2.4)	36
Table D-3. Config. Switches 6 - 10 (Ver. 2.4+)	36
Table D-4. USB-VCP interface serial port settings.	37
Table D-5. Microns/microsteps conversion.	38
Table D-6. Ranges.	38
Table D-7. TRIO MP-245 external control commands.	38

1. INTRODUCTION

1.1 Structure of the TRIO MP-245 Documentation Package

The TRIO MP-245 3-Axis Micromanipulator System is comprised of a ROE/controller, a power adapter, and a TRIO MP-245/M stepper-motor-based electromechanical micromanipulator. This manual consists of four parts: This chapter, Introduction, which provides an overview and general description of the TRIO MP-245 system; Chapter 2, Installation, which describes how to install, set up, and configure all components of the system; Chapter 3, Operations, which describes how to operate the TRIO MP-245; Chapter 4, Maintenance, describes how to perform routine and other maintenance; and Chapter 5, Reconfiguration, describes the reconfiguration possibilities of the TRIO MP-245 system.

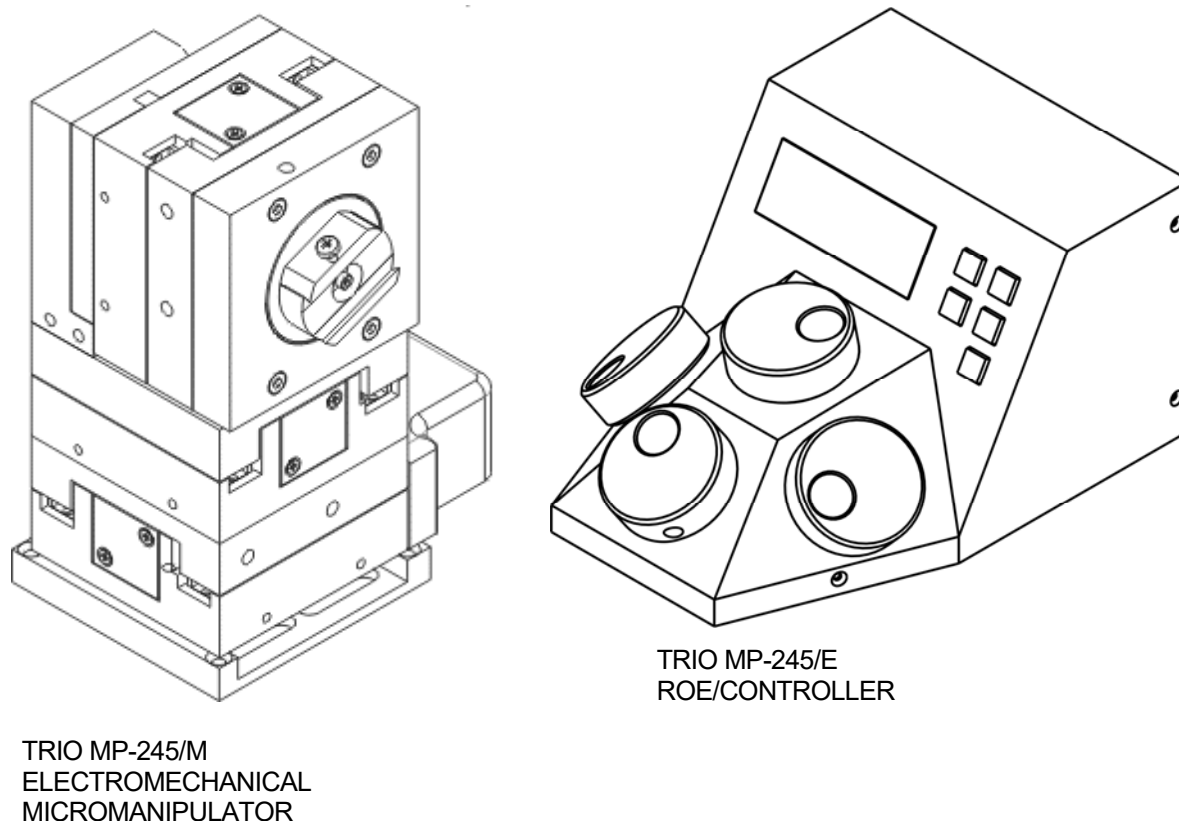


Figure 1-1. The TRIO MP-245 system

1.2 Components of the TRIO MP-245 System

Carefully remove all components from the shipping container. In addition to this manual, the following should be included:

- TRIO MP-245 ROE Rotary Optical Encoder input device with built-in controller and external power adapter.
- TRIO MP-245/M electromechanical micromanipulator
- 26-pin HD DSUB cable (connects the ROE/controller to the TRIO MP-245/M electromechanical micromanipulator).
- Power adapter
- Power adapter AC mains cable appropriate for your location
- Ground/Earth cable
- USB Cable

IMPORTANT

Once the TRIO MP-245 system has been unpacked, remove the shipping tape from the various locations on the TRIO MP-245/M electromechanical micromanipulator. The shipping tape must be removed before operating the TRIO MP-245 system. In the event that you need to transport the TRIO MP-245/M in the future, reapply 2 to 3-inch pieces of masking tape to the same locations. Once the tape has been removed, handle the TRIO MP-245/M with care. The mechanisms can be damaged if any of the axes are inadvertently moved without the tape in place.

1.3 Overview

1.3.1 Features

- Three independent axes (X, Y, and Z) each with 25mm travel with a virtual fourth axis (D) for coaxial pipette movement utilizing a tangent function factoring the holder's angle and the X and Z axes.
- Sub-micron 100nm resolution
- Digital display indicates coordinates in relative or absolute
- User-friendly, fanless compact controller with ROE preserves bench space
- Push button control of multiple functions – work, home, Lock, pulse and relative
- Robotic home and work position moves for easy automated pipette exchange

1.3.2 Description

The **TRIO MP-245**, the newest Sutter Instrument motorized manipulator, is easy to use and has three independent axes. The X, Y, and Z axes provide 25mm range of motion. D-axis movement is accomplished virtually using a tangent function of the chosen angle of the holder and simultaneously moving X and Z. The ROE controller has a digital display and keys for Home, Work, Pulse, Lock, and Relative. The compact, intuitive controller takes up minimal bench space, is fan-free, and easy to use.

While the axes provide X and Y orthogonal motion typical of most motorized manipulators, Sutter has introduced a diagonal axis with the **TRIO MP-245** so one can move the electrode coaxially at the exact desired angle of approach.

The **TRIO MP-245's** ROE provides fine control of electrode position and the rate of rotation of ROE dials for each axis determines the speed of travel. The finest step size is less than 100nm. Five conveniently located buttons on the ROE provide control of all the basic functions you will need in normal operation (Work, Home, Lock, Relative, and Pulse).

Press and hold **WORK** (for 3 seconds) to quickly store a work position, tap **HOME** to move all axes to an initial location that is useful for changing electrodes, or press and hold the **HOME** button (for 3 seconds) to memorize a new **HOME** position.

When ready to record data, the motor drive electronics can be suppressed by pressing the **LOCK** button. In the **LOCK** mode, the display turns red and ROE input is locked out to avoid any accidental motion.

Pressing and holding the **RELATIVE** button for three seconds at any location causes the display coordinates to all zeroes. When activating relative mode, the display turns blue.

To return to viewing the absolute coordinates, tap the **RELATIVE** button to toggle back. Finally, tapping the **PULSE** button causes a 3 μ m advance in the diagonal. This rapid burst of forward motion can assist in sharp electrode cell penetration.

All the electronics, except for a small power supply, are housed within the **TRIO MP-245** ROE and no separate controller or computer is required.

External computer control of the **TRIO MP-245** is possible via the USB connector mounted on the controller/ROE's rear panel. The controller's internal software is programmed with a defined set of commands allowing for a wide range of micromanipulator/stage movements as programmed in software residing in an external computer connected via USB.

2. INSTALLATION

When installing the **TRIO MP-245** system for the first time, it is recommended that the components of the system be installed in the following order: **TRIO MP-245/M** electromechanical micromanipulator first, followed by the **TRIO MP-245/E** ROE/Controller.

2.1 Mounting Instructions

The following sections describe how to mount the **TRIO MP-245/M** manipulator to a stand using the mounting adapter plate, how to adjust the pipette angle and how to mount different headstages.

2.1.1 Mounting the **TRIO MP-245/M** to the Stand or Platform

The **TRIO MP-245/M** attaches to the mounting adapter plate using four M3.5x6 hex head locking screws.

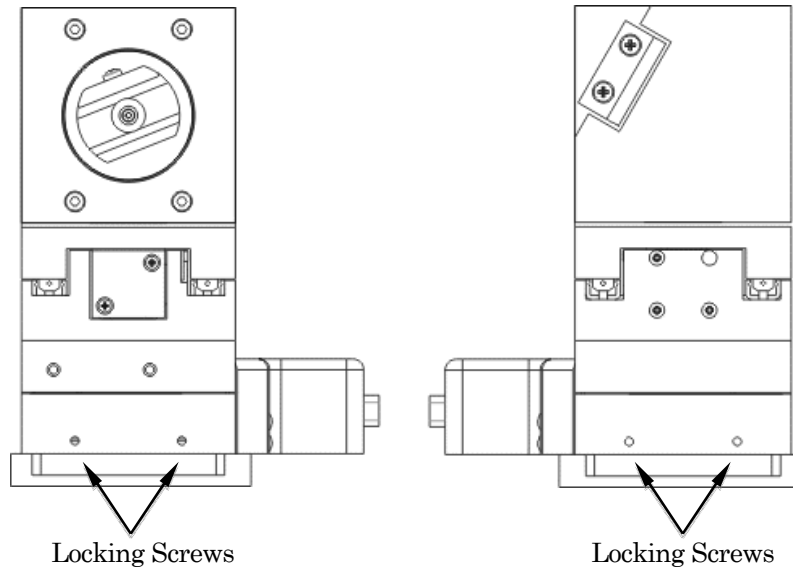


Figure 2-1. Side view of TRIO MP-245/M showing mounting adapter plate and lock screws.

The TRIO MP-245/M is shipped with the adapter plate in place. It is attached using four tapered pegs, along with four locking screws.

To remove it, first loosen the four hex screws that secure the manipulator to the pegs in the adapter plate. The rear pair is in a similar location in the back of the manipulator. Once the locking screws are sufficiently loosened, lift the TRIO MP-245/M upwards from the adapter plate.

Before attaching the adapter plate to the TRIO MP-245/M, you need to decide where to position the manipulator on your stand/platform. The stand can be any flat surface carrying $\frac{1}{4}$ -20, 10-32, or M6 holes on one-inch centers (such as a Sutter Instrument MT-series stand or MD series platform).

Examine the space of the platform onto which installation is to take place. Attach the control cable to TRIO MP-245/M and move the entire unit around on the platform until the precise desired position is determined. A small bag containing the necessary hardware to attach the TRIO MP-245 to the stand is included.

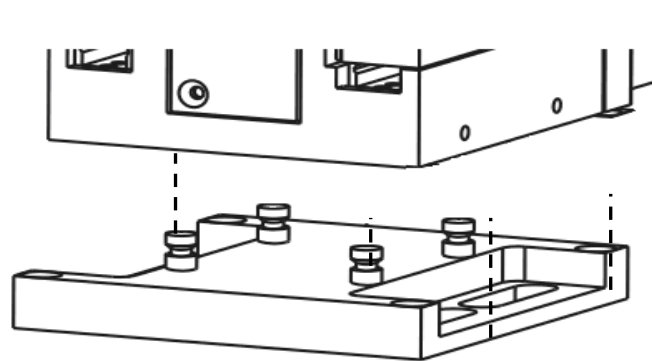


Figure 2-2. Mounting the TRIO MP-245/M on the Adapter Plate

Once the plate is mounted, align the pegs on top of the plate with the holes in the manipulator, push the X-axis firmly onto the plate, and re-tighten the locking hex set screws.

2.2 Headstage Mounting

Sutter IPA headstage, Axon headstages 203B or CV-7, and the Heka EPC-10 headstage have an integral dovetail that fits directly into the rotary dovetail slide bracket on the TRIO MP-245/M. The dovetail slide bracket on the TRIO MP-245/M also supports older Axon and Heka headstages when using the 4" dovetail extension.


Rod-mounted headstages and micro tools are accommodated by the use of a rod clamp that fits into the dovetail (not shown). All the headstage adapters and mounting hardware are included with the manipulator and are shipped in a zip lock plastic bag.

2.3 Other Accessories

One or more accessories may have been ordered and received for mounting the TRIO MP-245/M and/or modifying the headstage mount to the manipulator (i.e., rotating base, microscope stage mount, gantry, dovetail extension, etc.). Setup of these accessories is normally covered in documentation accompanying the accessory.

2.4 Electrical Connections and Initial Operating Instructions

Initially, you may want to simply connect the TRIO MP-245/M micromanipulator and the ROE/Controller together and try some gross movements in order to get a feel for the controls and how to make simple movements. It is perfectly acceptable to set the manipulators in the middle of a bench top, make all electrical connections and then observe each unit's movement by eye.

 **CAUTION:** Unless the TRIO MP-245/M micromanipulator electromechanical baseplate is firmly bolted down to a breadboard or solidly to a firm surface, the TRIO MP-245/M is likely to tip over when fully extending all of its axes, especially if its loaded with a headstage that extends beyond the TRIO MP-245/M's current center of gravity.

Upon deciding to directly install the TRIO MP-245 system in your rig, it is useful to follow the initial setup procedure to learn how to move the units to allow easy access to the mounting screws.

1. With the power switch on the back of the ROE in the OFF (0) position, connect the power adapter's 24VDC cable to the POWER receptacle.

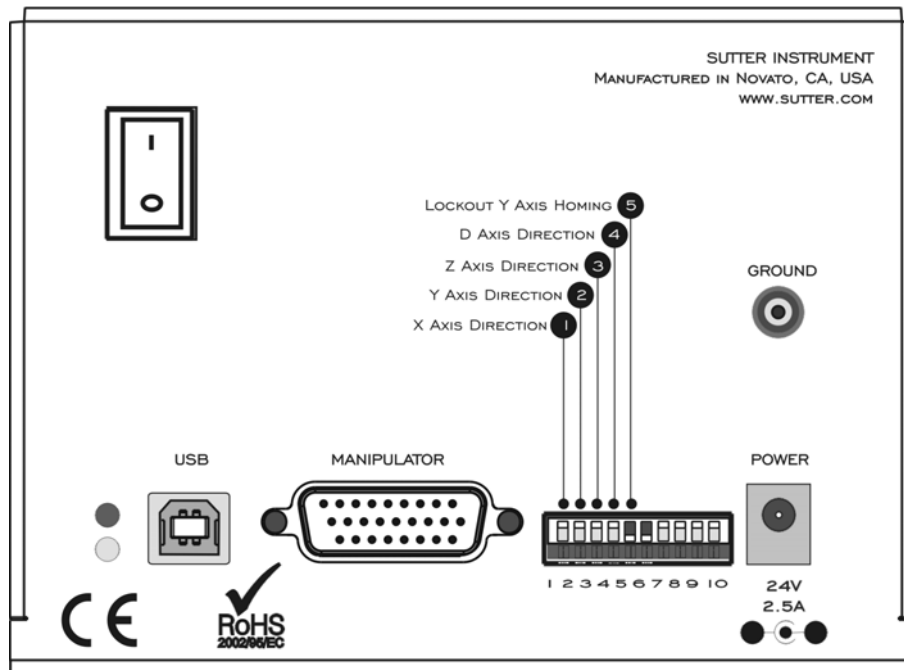


Figure 2-3. Rear of TRIO MP-245 ROE/Controller cabinet

2. With the power OFF (rear panel switch in the “0” position), connect a well-grounded/earthed wire to the GROUND banana plug receptacle.
3. With the power OFF, connect the male end of the DB-HD-26 cable to the MANIPULATOR connector on the ROE, the other end of which is connected to the TRIO MP-245/M micromanipulator electromechanical. (**See cautionary note below.**)
4. Verify that the six switches on the rear of the ROE are set as desired.
5. Power up the system by moving the power switch on the rear of the ROE to the “1” position.

*  **CAUTION: NEVER CONNECT OR DISCONNECT THE ROE/CONTROLLER FROM THE TRIO MP-245/M WHILE THE POWER IS ON!**

2.5 ROE/Controller Rear Panel Controls and Configuration

2.5.1 Power Switch

The power switch for the TRIO MP-245 system is located on the rear panel of the ROE/controller. At power up, the microprocessor in the ROE/controller scans the attached equipment and configures the system accordingly.

2.5.2 Configuration Switches

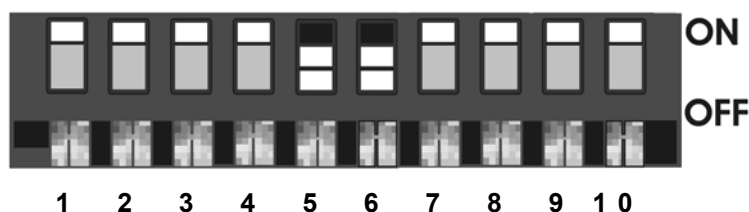


Figure 2-4. Configuration switches on rear of TRIO MP-245 ROE/Controller unit (switch positions shown are factory defaults).

2.5.2.1 Switches 1, 2, 3 and 4

These switches set the directionality for each of the four axes.

Table 2-1. Configuration Switches 1 - 4: Configuring the direction of each axis.

Switch #	Axis	Knob Rotation Directionality for Forward (+) Movement	
		Clockwise	Counterclockwise **
1	X	OFF (Up)*	ON (Down)
2	Y	OFF (Up)*	ON (Down)
3	Z	OFF (Up)*	ON (Down)
4	D	OFF (Up)*	ON (Down)

* Factory default (typical setting for right-hand-mounted manipulator).

** Possible (most often desired) setting for left-handed manipulator.

2.5.2.2 Switch 5: Y-Axis Lockout during Homing

Configures whether the Y axis is locked out while homing.

Table 2-2. Configuring the Homing Y-Movement Lock Out.

Switch #	Homing Y Movement Lock Out	
	Enabled	Disabled
5	OFF (Up)	ON (Down)*

* Factory default (recommended normal operation setting).

2.5.2.3 Switch 6: Sensor Test (Firmware < v2.2)

Table 2-3. Configuring the Sensor Test (Firmware < v2.2).

Switch #	Sensor Test	
	Enabled **	Disabled
6	OFF (Up)	ON (Down)*

* Factory default (do not change unless requested to by Sutter Instrument Technical Support).

**** CAUTION: To avoid damage to the micromanipulator or stage, DIP Switch 6 (Sensor Test) must always be set to ON (DOWN).**

2.5.2.4 Switch 6: Sensor Test (Firmware < v2.4)

Table 2-4. Configuration Switch 6: Configuring the Sensor Test.

Switch #	Calibration Homing on Power On	
	None (No calibration)	Calibrates to 1,000 μm for all axes at power on
6	OFF (Up)	ON (Down)*

*** CAUTION: Factory/normal default (do not change unless requested to by Sutter Instrument Technical Support).**

2.5.2.5 Switch 6: Calibration Homing on Power On (Firmware v2.4+)

Configures whether calibration homing occurs or not on power on (FW v2.4+).

Table 2-5. Configuration Switch 6: Calibration Homing on Power On (Firmware v2.4+).

Switch #	Calibration Homing on Power On	
	None (No calibration)	Calibrates to 1,000 μm for all axes at power on
6	OFF (Up)	ON (Down)*

* Factory default (recommended normal operation setting).

2.5.2.6 Switches 7 through 10

These four switches are reserved for all controllers programmed with a firmware version below v2.4. For firmware v2.4 and above, the switches are defined in the following paragraphs.

2.5.2.7 Switch 7: PULSE Button Functionality (Firmware v2.4+)

.

Table 2-6. Configuration Switch 7: PULSE button functionality.

Switch #	Definition	State	Setting	Position
7	PULSE button Mode: Pulse vs. Speed Select (0-3)	Speed-Select mode	OFF	UP
		Pulse mode	ON*	DOWN*

* Factory default (recommended normal operation setting)

2.5.2.8 Switch 8: Electromechanical Device Whole Step Size (Firmware v2.4+)

Table 2-7. Configuration Switch 8: Electromechanical device whole step size (1.5 vs. 2 μ m).

Switch #	Definition	State	Setting	Position
8	Electromechanical device whole step size	1.5 μ m Stepping (TRIO MP-245/M series)	OFF*	UP*
		Reserved	ON	DOWN

* Factory default (recommended normal operation setting)

2.5.2.9 Switch 9: Y-Axis Travel Length (Firmware v2.4+)

(FW v2.4)v2.4) This switch informs the ROE/controller as to the travel length of the Y axis in the connected micromanipulator electromechanical or stage. The length of the Y axis on the TRIO MP-245/M electromechanical micromanipulator is 25mm, so Switch 9 should be set to OFF (up) (factory default for a standard TRIO MP-245 system). For a variant model with 12.5mm of travel in the Y axis, Switch 9 must be set ON (down).

Table 2-8. Configuration Switch 9: Y-Axis travel length (Firmware v2.4+).

Switch #	Y-Axis Travel Length	
	25mm	12.5mm **
9	OFF (Up)*	ON (Down)

* Factory default (recommended normal operation setting).

** Use the ON (down) setting for half-length Y axis (12.5mm).

CAUTION: Always be certain that the position of Switch 9 correctly matches the physical length of travel of the Y axis on the connected device. Setting the switch to the OFF (up) position for 25mm could result in equipment damage if the attached device's Y axis length of travel is actually less than 25mm.

2.5.2.10 Switch 10: Length of X Axis (Firmware v2.4+)

(firmware v2.4)v2.4) This switch informs the ROE/controller as to the travel length of the X axis in the connected micromanipulator electromechanical or stage. The length of the X axis on the TRIO MP-245/M electromechanical micromanipulator is 25mm, so Switch 10 should

be set to OFF (up) (factory default for a standard TRIO MP-245 system). For a variant model with 50mm of travel in the X axis, Switch 10 must be set ON (down).

Table 2-9. Configuration Switch 10: Length of X Axis (Firmware v2.4+).

Switch #	Length of X Axis	
	25mm	50mm **
10	OFF (Up)*	ON (Down)

* Factory default (recommended normal operation setting).

** Use the ON (down) setting for double-length X axis (50mm).

CAUTION: Always be certain that the position of Switch 10 correctly matches the physical length of travel of the X axis on the connected device. Setting the switch to the ON (down) position for 50mm could result in equipment damage if the attached device's X axis length of travel is actually less than 50mm.

3. OPERATIONS

3.1 Main Controls and Indicators on the ROE/Controller


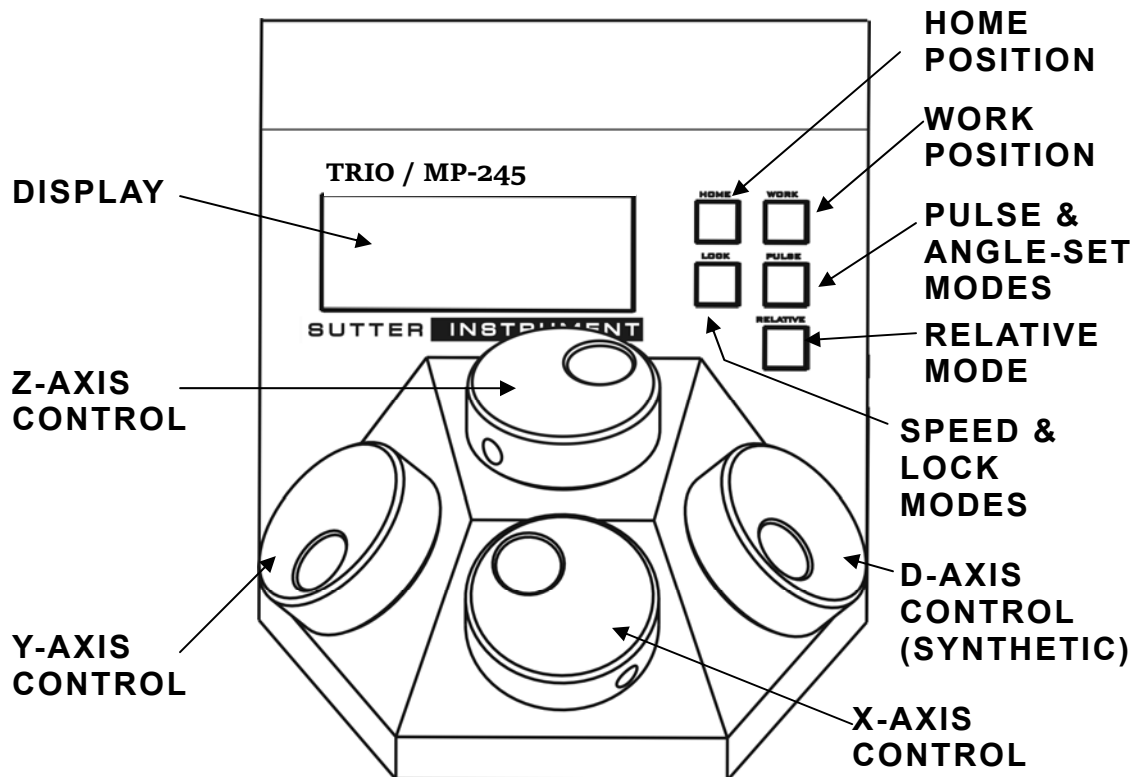


Figure 3-1. Front view of the TRIO MP-245 ROE/Controller

3.2 Display

3.2.1 Initial Startup

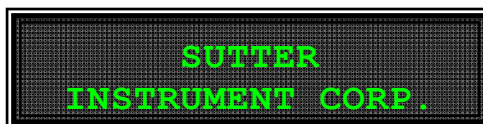

(Text in Green)


Figure 3-1. LCD Display showing startup screen.

When starting the TRIO MP-245 system for the first time or if the HOME position has not yet been defined (saved), the values of all four axes will be 1,000 micrometers (microns).



(Text in Green)

Figure 3-2. Startup screen

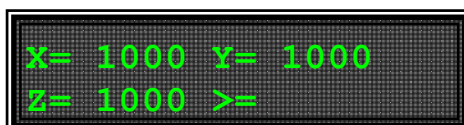

(Text in Green)

Figure 3-3. Factory default startup (Home) position


3.3 Control Operations

3.3.1 Maximum Positive Position Values:

Move the dial of an axis clockwise until its position value stops incrementing. The following table lists the maximum position value (in microns) for each axis.

Table 3-1. Maximum positive position value of each axis

Axis	Maximum Position Value (in microns)
X	25,000
Y	25,000
Z	25,000

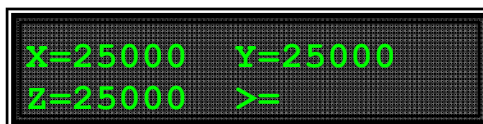

(Text in Green)

Figure 3-4. Maximum positive values

3.3.2 Setting Position for HOME or WORK

To set position, hold down HOME or WORK button for 3 seconds until beep sounds.

3.3.3 Setting the Angle of the Pipette/Headstage Holder

To change the angle of the holder, first loosen the set screw at the top of the rotary dovetail bracket, rotate the holder to the desired angle, and then retighten the set screw.

Measure the angle of the holder. (Tip: Many smart phones have an app with a level that can assist the user.)

3.3.4 Operating the Virtual D Axis

The TRIO MP-245 consists of three physical axes, X, Y and Z. A tangent function utilizing X and Z axes and the angle of the holder has been implemented to create a virtual D axis. Use between 10° and 90° for best results.

Zero (0°) is set with the diagonal being parallel to the table and 90° is set with the diagonal being perpendicular to the table.

To set the angle measured above, on the ROE hold LOCK down for several seconds. The screen will be red until the display indicates in green: “Select the angle in use (0-90)”.

Use the D dial on the ROE to set the value of the angle. Once this value is dialed in, do not touch the ROE knob for 8-10 sec. The virtual D angle will now be set.

3.3.5 Moving to the Home Position

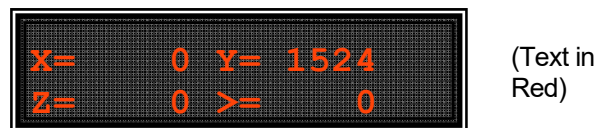


Figure 3-5. Moving to Home position (screen is amber while moving)

If the Home position has not yet been defined and saved, the Home position values for all axes will default to 1,000 microns, as shown in the following figure.

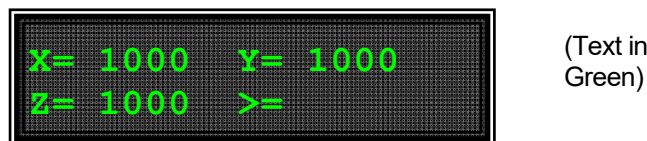


Figure 3-6. Factory default Home position

If the Home position has been previously defined (saved), pressing HOME will make a move to the defined home position (see example in the following figure).

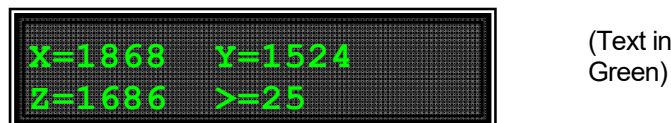


Figure 3-7. Example Home position defined and saved


To move to the Home position, press HOME. If the current position before pressing HOME is greater than the Home position, the movement will be as follows:

NOTE: Movement to the Home position works only if X coordinates of the HOME position are less than the WORK position.

1. Movement begins by retracting the Z axis (at the angle currently set) away from the sample.
2. Movement then continues along the X axis toward the Home position.
3. The final movement is along the Y-axis towards the operator and away from the microscope.

NOTE: Step 3 occurs only if “Y-Lockout” is disabled. Otherwise, no movement along the Y-axis occurs.

3.3.6 Moving to the Work Position

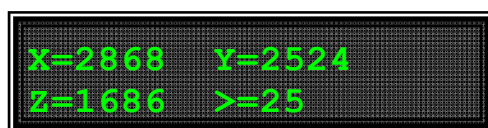

(Text in Green)

Figure 3-8. Example Work position

To move to the Work position, press the WORK button. If the current position before pressing WORK is less than the Work position, the movement will be as follows:


1. Movement travels along the Y-axis away from the operator and towards the microscope.
2. Movement is then made along the X axis toward the sample. Travel then continues along the diagonal until reaching its end-of-travel point.

NOTE: Step 1 occurs only if “Y-Lockout” is disabled. Otherwise, movement begins with Step 2.

3.3.7 Setting Absolute/Relative Coordinates Mode

The RELATIVE button toggles between Relative and Absolute coordinate systems. The default coordinate system on power up is Absolute, with the coordinates on the screen shown in green. To switch to relative coordinates, press the RELATIVE button once. To reset the current position to all zeroes, depress the RELATIVE for 3 seconds or until a beep is heard, and then release the button. This resets the current position to all zeroes.


Press RELATIVE once (briefly for < 2 sec.)

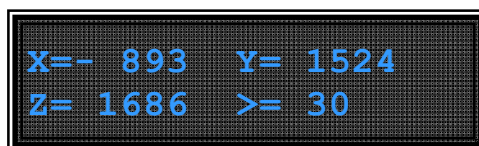

(Text in Blue)

Figure 3-9. Relative mode


Depress RELATIVE for 3 sec. or until beep sounds

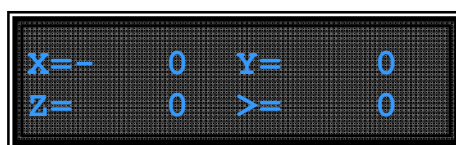

(Text in Blue)

Figure 3-10. Relative mode

Pressing RELATIVE briefly while in Relative mode, returns displayed coordinates back to Absolute mode

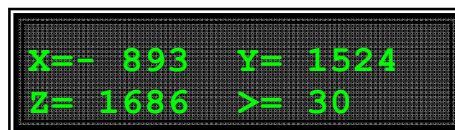

(Text in Green)

Figure 3-11. Absolute mode

3.3.8 Mode Indications

The TRIO MP-245 system has three modes of operation: Absolute coordinates, Relative coordinates, and Lock mode. The display turns color for each specific mode, as shown in the following table.

Table 3-2. Screen colors and modes

Screen Color	Mode	Example
Green	Absolute Coordinates	
Blue	Relative Coordinates	
Red	Knobs disabled during move to Home or Work position, while in Lock mode.	

3.3.9 Speed Control and ROE Knob Movements (SPEED)

The rate at which the ROE axis knobs move the electromechanical can be adjusted with the SPEED button. Each press of the button cycles through four speeds: 0 (normal) through 3 (fastest).

3.3.10 Movement Knobs Disabling and Lock Mode ([SPEED]/LOCK)

Axis-movement knobs are disabled during movements to Home, Work, or while in Lock Mode (display is in red).

3.3.11 Pausing Home Movements (HOME (while moving to Home))

After Move to Home has been initiated, and while the move is in progress, pressing HOME a second time pauses the manipulator. Pressing HOME again resumes movement.

3.3.12 Pausing Work Movements (WORK (while moving to Work))

After Move to Work has been initiated, and while the move is in progress, pressing WORK a second time pauses the manipulator. Pressing WORK again resumes movement.

3.3.13 Pulse Mode and Virtual D-Axis Movement (PULSE)

Pulse mode advances the D axis in $2.85\ \mu\text{m}$ steps. Each press of the PULSE button increments the Diagonal axis by one $2.85\text{-}\mu\text{m}$ step beyond the current position. This feature can be used to penetrate tough or resistant tissue.

3.4 Micropipette/Headstage Exchange

Mounted on the front of the Z-axis of the manipulator is the angle-control plate for the headstage mount.

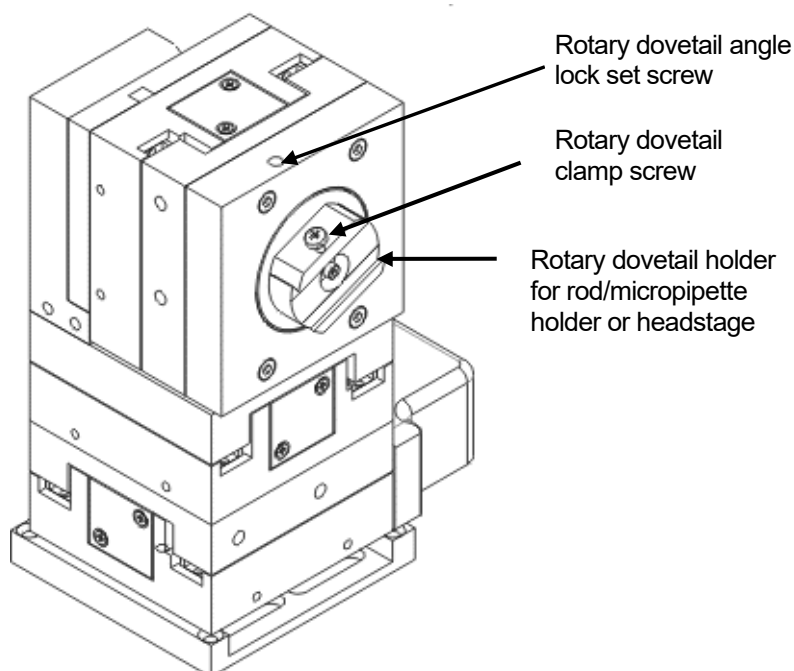


Figure 3-12. Angled side view of TRIO MP-245/M to change headstage mount

To change the headstage, loosen the screw in the center of the holding bracket. Slide the headstage upward out of the dovetail groove. Make any adjustments needed of the headstage, and then tighten down (but do not over tighten) the lock screw in the center of the holding bracket.

4. EXTERNAL CONTROL

4.1 General

Controlling the TRIO MP-245 externally via computer is accomplished by sending commands over the USB interface between the computer and the USB connector on the rear panel of the TRIO MP-245 controller/ROE. The USB device driver for Windows is downloadable from Sutter Instrument's web site (www.sutter.com). The TRIO MP-245 requires USB CDM (Combined Driver Model) Version 2.10.00 or higher. The CDM device driver for the TRIO MP-245 consists of two device drivers: 1) USB device driver, and 2) VCP (Virtual COM Port) device driver. Install the USB device driver first, followed by the VCP device driver. The VCP device driver provides a serial RS-232 I/O interface between a Windows application and the TRIO MP-245. Although the VCP device driver is optional, its

installation is recommended even if it is not going to be used. Once installed, the VCP can be enabled or disabled.

The CDM device driver package provides two I/O methodologies over which communications with the controller over USB can be conducted: 1) USB Direct (D2XX mode), or 2) Serial RS-232 asynchronous via the VCP device driver (VCP mode). The first method requires that the VCP device driver not be installed, or if installed, that it be disabled. The second method requires that the VCP be installed and enabled.

4.2 Virtual COM Port (VCP) Serial Port Settings

The following table lists the required RS-232 serial settings for the COM port (COM3, COM5, etc.) generated by the installation or enabling of the VCP device driver.

Table 4-1. USB-VCP interface serial port settings.

Property	Setting
Data (“Baud”) Rate (bits per second (bps))	57600
Data Bits	8
Stop Bits	1
Parity	None
Flow Control	None

The settings shown in the above table can be set in the device driver’s properties (via the Device Manager if in Windows) and/or programmatically in your application.

4.3 Protocol and Handshaking

Command sequences do not have terminators. All commands return an ASCII CR (Carriage Return; 13 decimal, 0D hexadecimal) to indicate that the task associated with the command has completed. When the controller completes the task associated with a command, it sends ASCII CR back to the host computer indicating that it is ready to receive a new command. If a command returns data, the last byte returned is the task-completed indicator.

4.4 Command Sequence Formatting

Each command sequence consists of at least one byte, the first of which is the “command byte”. Those commands that have parameters or arguments require a sequence of bytes that follow the command byte. No delimiters are used between command sequence arguments, and command sequence terminators are not used. Although most command bytes can be expressed as ASCII displayable/printable characters, the rest of a command sequence must generally be expressed as a sequence of unsigned byte values (0-255 decimal; 00 – FF hexadecimal, or 00000000 – 11111111 binary). Each byte in a command sequence being transmitted to the controller must contain an unsigned binary value. Attempting to code command sequences as “strings” is not advisable. Any command data being returned from the controller must also be received and initially treated as a sequence of unsigned byte values. Groups of contiguous bytes can later be combined to form larger values, as appropriate (e.g., 2 bytes into 16-bit “word”, or 4 bytes into a 32-bit “long” or “double word”). For the TRIO MP-245, all axis position values (number of microsteps) are stored as “unsigned long” (32-bit) values, and each is transmitted and received to and from the controller as four contiguous bytes.

4.5 Axis Position Command Parameters

All axis positional information is exchanged between the controller and the host computer in terms of microsteps. Conversion between microsteps and microns (micrometers) is the responsibility of the software running on the host computer. The number of microsteps for any axis position must always be exchanged between the controller and the application running on an external computer as an unsigned 32-bit value (“unsigned long” for C/C++ or “U32” for LabVIEW). “Unsigned” means the value is always positive; negative values are not allowed. An unsigned long or U32 consists of four contiguous bytes, with a byte/bit-ordering format of Little Endian (“Intel”) (most significant byte (MSB) in the first byte and least significant (LSB) in the last byte). If the platform on which your application is running is Little Endian, then no byte order reversal of axis position values is necessary. Examples of platforms using Little Endian formatting include any system using an Intel processor (including Microsoft Windows and Apple Mac OS X).

If the platform on which your application is running is Big Endian (e.g., Motorola PowerPC CPU), then these 32-bit position values must have their bytes reverse-ordered after receiving from, or before sending to, the controller. Examples of Big Endian platforms include many non-Intel-based systems, LabVIEW (regardless of operating system & CPU), and Java (programming language/environment). MATLAB adapts to the system on which it is running, so Little Endian may need to be enforced if running on a Big Endian system.

4.6 Microsteps and Microns (Micrometers)

All coordinates sent to and received from the controller are in microsteps. To convert between microsteps and microns (micrometers), use the following conversion factors (multipliers):

Table 4-2. Microns/microsteps conversion.

System/Device	From/To Units	Conversion Factor (multiplier)
TRIO MP-245 with TRIO MP-245/M micromanipulator	$\mu\text{steps} \rightarrow \mu\text{m}$	0.09375
	$\mu\text{m} \rightarrow \mu\text{steps}$	10.666666666667

For accuracy in your application, these conversion factors should be typed as double precision (“double”); “float” is not recommended. If the result is in microsteps, it can be typed as a 32-bit “long” integer; otherwise, it should be typed as floating point, preferably as double precision (“double”).

Table 4-3. Ranges (travel distances and bounds of each axis for compatible devices).

Device	Axis	Millimeters	Microns	Microsteps
TRIO MP-245/M	X, Y, & Z	0 – 25	0 – 25,000	0 – 266,667
With long X	X	0 – 50	0 – 50,000	0 – 533,334
With short Y	Y	0 – 12.5	0 – 12,500	0 – 133,334

4.7 Commands

4.7.1 Get Current Position (‘c’ or ‘C’) Command

This command is used to obtain the current position (X, Y, & Z coordinates) of the manipulator or stage and the current angle setting. The command sequence consists of one

byte as shown in the following table, followed by fourteen bytes containing X, Y, & Z position values in microsteps (4 bytes each), angle in degrees (1 byte), and completion indicator (1 byte).

Table 4-4. Get Current Position ('c' or 'C') command.

Command	Tx/ Delay/ Rx	Ver .	Total Bytes	Byte Offset (Len.)	Value			Alt- key- pad #	Ctrl- char	ASCII def./- char.	Description	
					Dec.	Hex.	Binary					
Get Current Position and Angle ('c' or 'C')	Tx	All	1	0	99	63	0110 0011	0099		'c' or 'C'	Returns the current positions (μ steps) of X, Y, & Z axes and angle setting (degrees).	
					or 67	or 43	or 0100 0011	or 0043				
	Rx	All	14	Three 4-byte (32-bit) values (current positions in μ steps of X, Y, & Z), + 1 byte for angle, + 1 byte for completion indicator. See Ranges table for minimum and maximum values.								
				0 (4)	X pos. in μ steps							
				4 (4)	Y pos. in μ steps							
				8 (4)	Z pos. in μ steps							
12				Angle in degrees								
13	13	0D	0000 1101		^M	<CR>	Completion indicator					

4.7.2 Move to Controller-Defined HOME Position ('h') Command

Table 4-5. Move to controller-defined HOME position ('h') command.

Command	Tx/ Delay/ Rx	Ver .	Total Bytes	Byte Offset (Len.)	Value			Alt- key- pad #	Ctrl- char	ASCII def./- char.	Description
					Dec.	Hex.	Binary				
Move to HOME Position ('h')	Tx		1	0	104	68	0110 1000	0104		'h'	Moves to the position saved for the controller's HOME button.
	Rx		1	0	13	0D	0000 1101			<CR>	Completion indicator

4.7.3 Move to Controller-Defined WORK Position ('w') Command

Table 4-6. Move to controller -defined WORK position ('w') command.

Command	Tx/ Delay/ Rx	Ver .	Total Bytes	Byte Offset (Len.)	Value			Alt- key- pad #	Ctrl- char	ASCII def./- char.	Description
					Dec.	Hex.	Binary				
Move to WORK Position ('w')	Tx		1	0	119	77	0111 0111	0119		'w'	Moves to the position saved for the controller's WORK button.
	Rx		1	0	13	0D	0000 1101			<CR>	Completion indicator

4.7.4 Move to Specified “Home” Position (‘H’) Command

This command instructs the controller to move all three axes to the position specified, moving Z axis first, followed by X and Y together. The command sequence consists of thirteen bytes.

Table 4-7. Move to specified “Home” position (‘H’) command.

Command	Tx/ Delay/ Rx	Ver .	Total Bytes	Byte Offset (Len.)	Value			Alt- key- pad #	Ctrl- char	ASCII def./- char.	Description		
					Dec.	Hex.	Binary						
Move to Specified “Home” Position (‘H’)	Tx	All	13	0	72	48	0100 1000	0072		‘H’	Move all 3 axes to specified position, moving Z before X & Y together (see Ranges table)		
				1 (4)	X μ steps								
				5 (4)	Y μ steps								
				9 (4)	Z μ steps								
	Rx	All	1	0	13	0D	0000 1101		^M	<CR>	Completion indicator		

4.7.5 Move to Specified “Work” Position (‘W’) Command

This command instructs the controller to move all three axes to the position specified, moving X and Y together, and then Z). The command sequence consists of thirteen bytes.

Table 4-8. Move to specified “Work” position (‘W’) command.

Command	Tx/ Delay/ Rx	Ver .	Total Bytes	Byte Offset (Len.)	Value			Alt- key- pad #	Ctrl- char	ASCII def./- char.	Description		
					Dec.	Hex.	Binary						
Move to Specified “Work” Position (‘W’)	Tx	All	13	0	87	57	0101 0111	0087		‘W’	Move all 3 axes to specified position, moving X & Y together first, and then Z (see Ranges table)		
				1 (4)	X μ steps								
				5 (4)	Y μ steps								
				9 (4)	Z μ steps								
	Rx	All	1	0	13	0D	0000 1101		^M	<CR>	Completion indicator		

4.7.6 Move in Straight Line to Specified Position at Specified Speed ('S') Command

This command instructs the controller to move all three axes in a straight line to the position specified (X, Y, & Z simultaneously). The command sequence consists of seventeen bytes.

Table 4-9. Straight-line move to specified position ('S') command.

Command	Tx/ Delay/ Rx	Ver .	Total Bytes	Byte Offset (Len.)	Value			Alt- key- pad #	Ctrl- char	ASCII def./- char.	Description
					Dec.	Hex.	Binary				
Move in Straight Line to Specified Position at Specified Speed ('S')	Tx	All	14	0	83	53	1001 0111	0083		'S'	Move all three axes in a straight line to specified position (see Ranges table)
				1	15	0F	0000 1111	0015	^O	Speed (15 – 0 (fastest through slowest))	
					-	-	-	-	-		
				2 (4)						X μ steps	
				6 (4)						Y μ steps	
			10 (4)				Z μ steps				
	Rx	All	1	0	13	0D	0000 1101		^M	<CR>	Completion indicator

4.7.7 Interrupt Straight-Line Move ('^C') Command

This command interrupts a move in progress only if the move was initiated by the “move in straight line ('S') command. The command sequence consists of one byte.

Table 4-10. Interrupt a straight-line move in progress ('^C') command.

Command	Tx/ Delay/ Rx	Ver .	Total Bytes	Byte Offset (Len.)	Value			Alt- key- pad #	Ctrl- char	ASCII def./- char.	Description
					Dec.	Hex.	Binary				
Interrupt Straight-Line Move (^C)	Tx	All	1	0	3	03	0000 0011	0003	^C	<ETX>	Interrupts a move in progress (only for moves initiated by the “Straight-line” move ('S') command)
				Rx		1	0	13	0D	0000 1101	

4.7.8 Move to Specified X-Axis Position ('x' or 'X') Command

This command moves to a specified position for only the X-axis.

Table 4-11. Move to specified X-axis position ('x' or 'X') command.

Command	Tx/ Delay/ Rx	Ver .	Total Bytes	Byte Offset (Len.)	Value			Alt- key- pad #	Ctrl- char	ASCII def./- char.	Description
					Dec.	Hex.	Binary				
Move to specified X axis Position (‘x’ or ‘X’)	Tx	All	5	0	120	78	0111 1000	0120		'x'	Move X axis to specified position (see Ranges table)
					or	or	or	or		or	
				90	5A	0101 1010	0090		'X'		
			1 (4)							X μ steps	
	Rx		1	0	13	0D	0000 1101			<CR>	Completion indicator

4.7.9 Move to Specified Y-Axis Position ('y' or 'Y') Command

This command moves to a specified position for only the Y-axis.

Table 4-12. Move to specified Y-axis position ('y' or 'Y') command.

Command	Tx/ Delay/ Rx	Ver .	Total Bytes	Byte Offset (Len.)	Value			Alt- key- pad #	Ctrl- char	ASCII def./- char.	Description
					Dec.	Hex.	Binary				
Move to specified Y axis Position (‘y’ or ‘Y’)	Tx	All	5	0	121 or 91	79 or 5B	0111 1001 or 0101 1011	0121 or 0091		'y' or 'Y'	Move Y axis to specified position (see Ranges table)
					1 (4)					Y μ steps	
	Rx		1	0	13	0D	0000 1101			<CR>	Completion indicator

4.7.10 Move to Specified Z-Axis Position ('z' or 'Z') Command

This command moves to a specified position for only the Z-axis.

Table 4-13. Move to specified Z-axis position ('z') command.

Command	Tx/ Delay/ Rx	Ver .	Total Bytes	Byte Offset (Len.)	Value			Alt- key- pad #	Ctrl- char	ASCII def./- char.	Description
					Dec.	Hex.	Binary				
Move to specified Z axis Position (‘z’ or ‘Z’)	Tx	All	5	0	122 or 92	7A or 5C	0111 1010 or 0101 1100	0122 or 0092		'z' or 'Z'	Move Z-axis to specified position (see Ranges table)
					1 (4)					Z μ steps	
	Rx		1	0	13	0D	0000 1101			<CR>	Completion indicator

4.7.11 Setting the Angle ('A') Command

Sets the angle value, in degrees, to match the angle position of the rotary dovetail

Table 4-14. Set the angle ('A') command.

Command	Tx/ Delay/ Rx	Ver .	Total Bytes	Byte Offset (Len.)	Value			Alt- key- pad #	Ctrl- char	ASCII def./- char.	Description
					Dec.	Hex.	Binary				
Enter Angle (‘A’)	Tx	All	2	0	65	41	1010 1001	0065		'A'	Sets the angle value, in degrees, to match the angle position of the rotary dovetail
					1	0 - 90	00 - 5A	0000 0000 - 0101 1010	0000 - 0090	<NULL> - 'z'	Angle in degrees between 0 and 90.
	Rx		1	0	13	0D	0000 1101			<CR>	Completion indicator

4.7.12 Notes

1. All positions sent to and received from the controller are in microsteps. See Microns/microsteps conversion for conversion between microns (micrometers (μm)) and microsteps (μsteps).

2. See Ranges for exact minimum and maximum values for each axis of each compatible device that can be connected.
3. A short delay (usually around 2 ms) is recommended between commands (after the reception of one command and the sending of the next command).
4. All positions sent and received to and from the controller are in microsteps and consist of 32-bit positive integer values (four contiguous bytes). For C/C++, these are typed as “unsigned long”; “U32” for LabVIEW. Although a signed type might be desirable for computation and compatibility purposes, the application must ensure that only unsigned longs or U32 values are sent to, and received from, the controller (i.e., no negative values).
5. All 32-bit position values transmitted to, and received from, the controller must be bit/byte-ordered in “Little Endian” format. This means that the least significant bit/byte is last (last to send and last to receive). Byte-order reversal may be required on some platforms. Microsoft Windows, Intel-based Apple Macintosh systems running Mac OS X, and some Intel/AMD processor based Linux distributions handle byte storage in Little-Endian byte order so byte reordering is not necessary before converting to/from 32-bit “long” values. LabVIEW always handles “byte strings” in “Big Endian” byte order irrespective of operating system and CPU, requiring that the four bytes containing a microsteps value be reverse-ordered before/after conversion to/from a multibyte type value (I32, U32, etc.). MATLAB automatically adjusts the endianness of multibyte storage entities to that of the system on which it is running, so explicit byte reordering is generally unnecessary unless the underlying platform is Big Endian. If your development platform does not have built-in Little/Big Endian conversion functions, bit reordering can be accomplished by first swapping positions of the two bytes in each 16-bit half of the 32-bit value, and then swap positions of the two halves. This method efficiently and quickly changes the bit ordering of any multibyte value between the two Endian formats (if Big Endian, it becomes Little Endian, and if Little Endian, it becomes then Big Endian).
6. “Move” commands might have short to long distances of travel. If not polling for return data, an appropriate delay should be inserted between the sending of the command sequence and reception of return data so that the next command is sent only after the move is complete. This delay can be auto-calculated by determining the distance of travel (difference between current and target positions) and rate of travel. This delay is not needed if polling for return data. In either case, however, an appropriate timeout must be set for the reception of data so that the I/O does not time out before the move is made and/or the delay expires.
7. Actual speed for the “Straight-Line Move ‘S’” command can be determined with the following formula: $(5000 / 16) * (sp + 1)$, where 5000 is the maximum speed in microns/second and “sp” is the speed level 0 (slowest) through 15 (fastest). For mm/second or microns/millisecond, multiply result by 0.001.

5. MAINTENANCE

Routine cleaning of the TRIO MP-245 system is required to prevent excessive dust accumulations. Wipe all exterior surfaces with a dry, soft, cotton cloth.

Periodically inspect all cables and connections to make sure that all connections are made well and that all connectors are well and evenly seated.

6. RECONFIGURATION

6.1 Changing the Rotary Knob Functions on the ROE/Controller

The axis motor assignment of each axis control knob on the ROE can be changed by opening the ROE/Controller cabinet as seen in the figure below and changing cables to appropriate connectors.

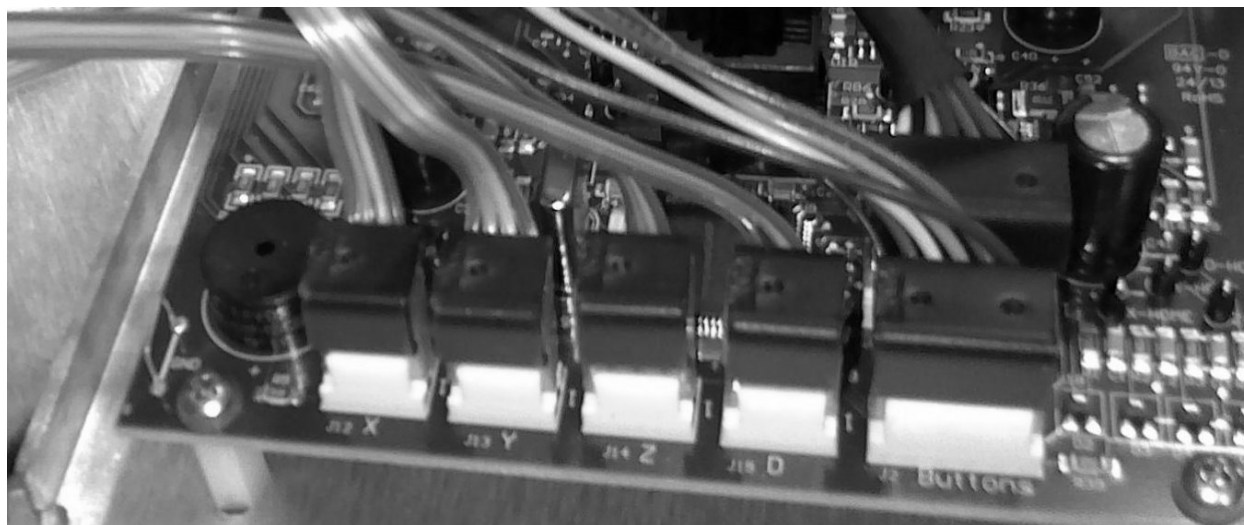


Figure 6-1. Locations of the axis connectors inside the ROE/Controller

APPENDIX A. LIMITED WARRANTY

- Sutter Instrument Company, a division of Sutter Instrument Corporation, limits the warranty on this instrument to repair and replacement of defective components for two years from date of shipment, provided the instrument has been operated in accordance with the instructions outlined in this manual.
- Abuse, misuse, or unauthorized repairs will void this warranty.
- Warranty work will be performed only at the factory.
- The cost of shipment both ways is paid for by Sutter Instrument during the first three months this warranty is in effect, after which the cost is the responsibility of the customer.
- The limited warranty is as stated above and no implied or inferred liability for direct or consequential damages is intended.
- An extended warranty for up to three additional years can be purchased at the time of ordering, or until the original warranty expires. For pricing and other information, please contact Sutter Instrument.

APPENDIX B. ACCESSORIES

W621 150	Ground cable
285204	4 inch dovetail extension
285210	Mounting adapter plate
225RBI	Rotating base
221165	Z-axis vertical extension
BR-AW	Rod holding clamp for XenoWorks ^{®*} injectors (for rod OD 2-4 mm)
MP-ROD	Rod holder (for rod OD 6.25 mm or larger)
MD Series	Micromanipulator platform
MT-78-FS	Large fixed-stage platform
MT-78-FS/M6	Large fixed-stage platform with M6 tapped holes
MT-75	Standard gantry stand 8.7 to 13.4 in (22.1 to 34 cm)
MT-75S	Short gantry stand 6.7 to 9.6 in (17 to 24.4 cm)
MT-75T	Tall gantry stand 10.7 to 15.4 in (27.2 to 39.1 cm)
MT-75XT	Extra tall gantry stand 14.7 to 18.5 in (37.4 to 47 cm)
MP-RISER-0.5[†]	½-inch (12.5mm) riser
MP-RISER-1.0	1-inch (25mm) riser

* XenoWorks[®] is a registered trademark of Sutter Instrument Company.

† Risers may be combined to achieve desired height

APPENDIX C. TECHNICAL SPECIFICATIONS




Travel	25mm on X, Y, and Z axes
Resolution	Minimal microstep size is 62.5 nanometers per microstep. Display has single micron resolution.
Speed	5 mm/sec. maximum
Drift	< 1 micron/hr drive mechanism
Electrical:	
Power Adapter:	Meanwell GS60A24-P1J
Input (mains)	100 - 240 VAC, 50/60 Hz, 1.4A
Output (to controller)	24V DC, 2.5A, 60W Max. (see following table for cable specs)
System Power consumption	60-Watts maximum
Mains fuses	None replaceable (power protection built into the Power Adapter)
Cables	(Refer to the following tables for a description of all possible cables.)

Table C-1. TRIO MP-245 cables and receptacles/connectors.

Controller Rear Panel Port Connector/-Receptacle	Cable Connector Types	Connects to ...	Cable Type	Cable Max. Length
Power Adapter 3-pin male connector	<p>←3-pin power standard (female)</p> <p>3-pin male→</p> <p>(Geographical region dependent)</p>	Mains power source.	10A, 250V, with safety ground plug	3 meters (approx. 10 feet)
ROE/Controller Cabinet: MANIPULATOR (26-Pin HD DSUB female receptacle)	<p>←HD DB-26 male</p> <p>HD DB-26 female→</p> <p>(Straight-through)</p>	TRIO MP-245/M	Minimum of 26 awg stranded wire with 500 Volt.	3 meters (approx. 10 feet)
Power Adapter	<p>←(fixed)</p> <p>ID 2.1 x OD 5.5 mm Barrel Plug (male)→</p>	ROE/Controller Cabinet: POWER receptacle (center pin positive)	UL1185 18AWG	1.8 meters (approx. 6 feet)
ROE/Controller Cabinet: GROUND (1-pin Banana-style female receptacle)	<p>←Banana male</p> <p>→ Alligator clip (hooded)</p>	a ground/-earth source (user determined)		
ROE/Controller Cabinet: USB	<p>←A</p> <p>→B</p>	Computer USB port		

APPENDIX D. QUICK REFERENCE

D.1. Manual Operation

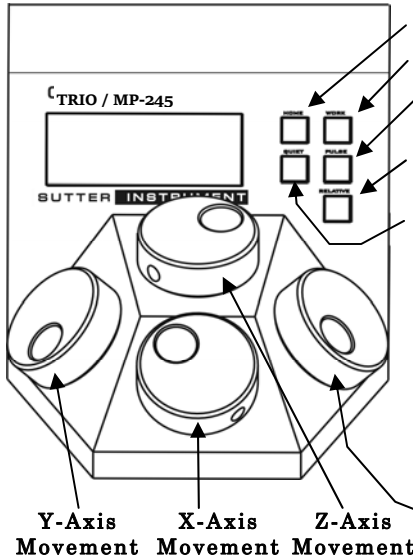

HOME: Move to defined home position. Press again to pause/resume.

WORK: Move to defined work position. Press again to pause/resume.

PULSE: Advances diagonal axis in 2.85 μm steps. Hold 3-sec. to enter **ANGLE SET** mode (active for 10 sec.) – use Knob D to change angle (0 – 90 degrees).

RELATIVE: Toggles between **Relative** and **Absolute** position moves. Hold 3-sec. to set the relative mode origin to the current absolute position.

SPEED: Cycles through Speed 0 (normal) through 3. Hold 3 sec. to enter **LOCK** mode.

Setting Home/Work Pos., Relative Mode Origin Pos., & Angle: To set position, hold down HOME, WORK, & RELATIVE buttons for 3 seconds until beep sounds. Likewise, set angle with LOCK.

Screen-color mode indications: Green = Absolute position; Blue = Relative position; Red = Movement in progress or in quiet (**LOCK**) mode; knobs disabled.

Movement Knobs Disabling and Quiet (LOCK**) Mode:** Movement knobs are disabled during movement to Home, Work, external movement command, or while in quiet (**LOCK**) Mode.

D-Axis Movement (Synthetic); select angle in degrees when in **ANGLE SET** mode.

D.2. Configuration

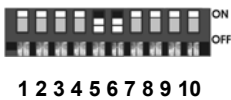
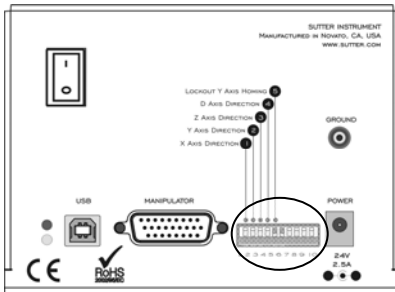


Table D-1. Configuration Switches 1 – 5.

Sw #	Definition	State	Setting	Position
1	X-Axis Knob Rotation for Forward (+) Movement	Clockwise	Off*	Up*
		Counter	On	Down
2	Y-Axis Knob Rotation for Forward (+) Movement	Clockwise	Off*	Up*
		Counter	On	Down
3	Z-Axis Knob Rotation for Forward (+) Movement	Clockwise	Off*	Up*
		Counter	On	Down
4	D-Axis Knob Rotation for Forward (+) Movement	Clockwise	Off*	Up*
		Counter	On	Down
5	Y Axis Lock Out for Homing	Enabled	Off	Up
		Disabled	On*	Down*

* Normal operation (factory default).

Table D-2. Configuration Switches 6 – 10 (Ver. <2.4)

Sw #	Definition	State	Setting	Position
6	Sensor Test** (see Caution)	Enabled	Off	Up
		Disabled	On*	Down*
7 - 10	Reserved		Off*	Up*

**CAUTION: To avoid damage to the micromanipulator or stage, DIP Switch 6 (Sensor Test) must always be set to OFF (DOWN).

Table D-3. Config. Switches 6 - 10 (Ver. 2.4+)

Sw #	Definition	State	Setting	Position
6	Calibration Homing on Power On	None	Off	Up
		Calibrate	On	Down*
7	SPEED Select or PULSE button mode	Speed Select	Off*	Up*
		Pulse	On	Down
8	Electromechanical device whole step size	1.5 μm	Off*	Up*
		Reserved	On	Down
9 ***	X-axis travel length	25mm	Off*	Up*
		50mm	On	Down
10 ***	Y-axis travel length:	25mm	Off*	Up*
		12.5mm	On	Down

***CAUTION: To avoid possible equipment damage, always be certain that Switch 9 and 10 settings correctly match the physical length of travel of the X & Y axes on the connected device.

D.3. External Control

Controlling the TRIO MP-245 externally via computer is accomplished by sending commands over the USB interface between the computer and the USB connector on the rear panel of the TRIO MP-245 controller/ROE. The USB device driver for Windows is downloadable from Sutter Instrument's web site (www.sutter.com). The TRIO MP-245 requires USB CDM (Combined Driver Model) Version 2.10.00 or higher. The CDM device driver for the TRIO MP-245 consists of two device drivers: 1) USB device driver, and 2) VCP (Virtual COM Port) device driver. Install the USB device driver first, followed by the VCP device driver. The VCP device driver provides a serial RS-232 I/O interface between a Windows application and the TRIO MP-245. Although the VCP device driver is optional, its installation is recommended even if it is not going to be used. Once installed, the VCP can be enabled or disabled.

The CDM device driver package provides two I/O methodologies over which communications with the controller over USB can be conducted: 1) USB Direct (D2XX mode), or 2) Serial RS-232 asynchronous via the VCP device driver (VCP mode). The first method requires that the VCP device driver not be installed, or if installed, that it be disabled. The second method requires that the VCP be installed and enabled.

Virtual COM Port (VCP) Serial Port Settings: The following table lists the required RS-232 serial settings for the COM port (COM3, COM5, etc.) generated by the installation or enabling of the VCP device driver.

Table D-4. USB-VCP interface serial port settings.

Property	Setting
Data ("Baud") Rate (bits per second (bps))	57600
Data Bits	8
Stop Bits	1
Parity	None
Flow Control	None

The settings shown in the above table can be set in the device driver's properties (via the Device Manager if in Windows) and/or programmatically in your application.

Protocol and Handshaking: Command sequences do not have terminators. All commands return an ASCII CR (Carriage Return; 13 decimal, 0D hexadecimal) to indicate that the task associated with the command has completed. When the controller completes the task associated with a command, it sends ASCII CR back to the host computer indicating that it is ready to receive a new

command. If a command returns data, the last byte returned is the task-completed indicator.

Command Sequence Formatting: Each command sequence consists of at least one byte, the first of which is the "command byte". Those commands that have parameters or arguments require a sequence of bytes that follow the command byte. No delimiters are used between command sequence arguments, and command sequence terminators are not used. Although most command bytes can be expressed as ASCII displayable/printable characters, the rest of a command sequence must generally be expressed as a sequence of unsigned byte values (0-255 decimal; 00 – FF hexadecimal, or 00000000 – 11111111 binary). Each byte in a command sequence being transmitted to the controller must contain an unsigned binary value. Attempting to code command sequences as "strings" is not advisable. Any command data being returned from the controller must also be received and initially treated as a sequence of unsigned byte values. Groups of contiguous bytes can later be combined to form larger values, as appropriate (e.g., 2 bytes into 16-bit "word", or 4 bytes into a 32-bit "long" or "double word"). For the TRIO MP-245, all axis position values (number of microsteps) are stored as "unsigned long" (32-bit) values, and each is transmitted and received to and from the controller as four contiguous bytes.

Axis Position Command Parameters: All axis positional information is exchanged between the controller and the host computer in terms of microsteps. Conversion between microsteps and microns (micrometers) is the responsibility of the software running on the host computer. The number of microsteps for any axis position must always be exchanged between the controller and the application running on an external computer as an unsigned 32-bit value ("unsigned long" for C/C++ or "U32" for LabVIEW). "Unsigned" means the value is always positive; negative values are not allowed. An unsigned long or U32 consists of four contiguous bytes, with a byte/bit-ordering format of Little Endian ("Intel") (most significant byte (MSB) in the first byte and least significant (LSB) in the last byte). If the platform on which your application is running is Little Endian, then no byte order reversal of axis position values is necessary. Examples of platforms using Little Endian formatting include any system using an Intel processor (including Microsoft Windows and Apple Mac OS X).

If the platform on which your application is running is Big Endian (e.g., Motorola PowerPC CPU), then these 32-bit position values must have their bytes reverse-ordered after receiving from, or before sending to, the controller. Examples of Big Endian

platforms include many non-Intel-based systems, LabVIEW (regardless of operating system & CPU), and Java (programming language/environment). MATLAB adapts to the system on which it is running, so Little Endian may need to be enforced if running on a Big Endian system.

Microsteps and Microns (Micrometers): All coordinates sent to and received from the controller are in microsteps. To convert between microsteps and microns (micrometers), use the following conversion factors (multipliers):

Table D-5. Microns/microsteps conversion.

System/Device	From/To Units	Conversion Factor (multiplier)
TRIO MP-245 with TRIO MP-245/M micromanipulator	$\mu\text{steps} \rightarrow \mu\text{m}$	0.09375
	$\mu\text{m} \rightarrow \mu\text{steps}$	10.666666666667

For accuracy in your application, these conversion factors should be typed as double precision

(“double”); “float” is not recommended. If the result is in microsteps, it can be typed as a 32-bit “long” integer; otherwise, it should be typed as floating point, preferably as double precision (“double”).

Table D-6. Ranges

Device	Axis	Length	Microns	Microsteps
TRIO MP-245/M	X, Y, & Z	25mm	0 – 25,000	0 – 266,667
With long X	X	50mm	0 – 50,000	0 – 533,334
With short Y	Y	12.5mm	0 – 12,500	0 – 133,334

Command Reference: The following table lists all the external-control commands for the TRIO MP-245.

Table D-7. TRIO MP-245 external control commands.

Command	Tx/-Delay/-Rx	Ver.	Total Bytes	Byte Offset (Len.)	Value			Alt-key-pad #	Ctrl-char	ASCII def./-char.	Description
					Dec.	Hex.	Binary				
Get Current Position and Angle ('c' or 'C')	Tx	All	1	0	99 or 67	63 or 43	0110 0011 or 0100 0011	0099 or 0043		'c' or 'C'	Returns the current positions (μsteps) of X, Y, & Z axes and angle setting (degrees).
	Rx	All	14		Three 4-byte (32-bit) values (current positions in μsteps of X, Y, & Z), + 1 byte for angle, + 1 byte for completion indicator. See Ranges table for minimum and maximum values.						
				0 (4)	X pos. in μsteps						
				4 (4)	Y pos. in μsteps						
				8 (4)	Z pos. in μsteps						
				12	Angle in degrees						
			13	13	0D	0000 1101		^M	<CR>	Completion indicator	
Move to HOME Position ('h')	Tx	All	1	0	104	68	0110 1000	0104		'h'	Moves to the position saved for the controller's HOME button.
	Rx	All	1	0	13	0D	0000 1101			<CR>	Completion indicator
Move to WORK Position ('w')	Tx	All	1	0	119	77	0111 0111	0119		'w'	Moves to the position saved for the controller's WORK button.
	Rx	All	1	0	13	0D	0000 1101			<CR>	Completion indicator
Move to Specified "Home" Position ('H')	Tx	All	13	0	72	48	0100 1000	0072		'H'	Move all 3 axes to specified position, moving Z before X & Y together (see Ranges table)
				1 (4)	X μsteps						
				5 (4)	Y μsteps						
				9 (4)	Z μsteps						
	Rx	All	1	0	13	0D	0000 1101		^M	<CR>	Completion indicator

Command	Tx/ Delay/ Rx	Ver .	Total Bytes	Byte Offset (Len.)	Value			Alt- key- pad #	Ctrl- char	ASCII def./- char.	Description	
					Dec.	Hex.	Binary					
Move to Specified "Work" Position ('W')	Tx	All	13	0	87	57	0101 0111	0087		'W'	Move all 3 axes to specified position, moving X & Y together first, and then Z (see Ranges table)	
				1 (4)	X μ steps							
				5 (4)	Y μ steps							
				9 (4)	Z μ steps							
	Rx	All	1	0	13	0D	0000 1101		^M	<CR>	Completion indicator	
Move in Straight Line to Specified Position at Specified Speed ('S')	Tx	All	14	0	83	53	1001 0111	0083		'S'	Move all three axes in a straight line to specified position (see Ranges table)	
				1	15 - 0	0F - 00	0000 1111 - 0000 0000	0015 - 0000	^O - ^@		Speed (15 – 0 (fastest through slowest))	
				2 (4)	X μ steps							
				6 (4)	Y μ steps							
				10 (4)	Z μ steps							
	Rx	All	1	0	13	0D	0000 1101		^M	<CR>	Completion indicator	
Interrupt Straight-Line Move (^C)	Tx	All	1	0	3	03	0000 0011	0003	^C	<ETX>	Interrupts a move in progress (only for moves initiated by the "Straight-line" move ('S') command)	
	Rx	All	1	0	13	0D	0000 1101			<CR>	Completion indicator	
Move to specified X axis Position ('x' or 'X')	Tx	All	5	0	120 or 90	78 or 5A	0111 1000 or 0101 1010	0120 or 0090		'x' or 'X'	Move X axis to specified position (see Ranges table)	
				1 (4)	X μ steps							
	Rx	All	1	0	13	0D	0000 1101			<CR>	Completion indicator	
Move to specified Y axis Position ('y' or 'Y')	Tx	All	5	0	121 or 91	79 or 5B	0111 1001 or 0101 1011	0121 or 0091		'y' or 'Y'	Move Y axis to specified position (see Ranges table)	
				1 (4)	Y μ steps							
	Rx	All	1	0	13	0D	0000 1101			<CR>	Completion indicator	
Move to specified Z axis Position ('z' or 'Z')	Tx	All	5	0	122 or 92	7A or 5C	0111 1010 or 0101 1100	0122 or 0092		'z' or 'Z'	Move Z-axis to specified position (see Ranges table)	
				1 (4)	Z μ steps							
	Rx	All	1	0	13	0D	0000 1101			<CR>	Completion indicator	
Enter Angle ('A')	Tx	All	2	0	65	41	1010 1001	0065		'A'	Sets the angle value, in degrees, to match the angle position of the rotary dovetail	
				1	0 - 90	00 - 5A	0000 0000 - 0101 1010	0000 - 0090		<NUL> - 'z'	Angle in degrees between 0 and 90.	
					Rx	All	1	0	13	0D	0000 1101	

NOTES:

1. All positions sent to and received from the controller are in microsteps. See **Microns/microsteps conversion** for conversion between microns (micrometers (μm)) and microsteps (μsteps).
2. See **Ranges** for exact minimum and maximum values for each axis of each compatible device that can be connected.
3. A short delay (usually around 2 ms) is recommended between commands (after the reception of one command and the sending of the next command).
4. All positions sent and received to and from the controller are in microsteps and consist of 32-bit positive integer values (four contiguous bytes). For C/C++, these are typed as “unsigned long”; “U32” for LabVIEW. Although a signed type might be desirable for computation and compatibility purposes, the application must ensure that only unsigned longs or U32 values are sent to, and received from, the controller (i.e., no negative values).
5. All 32-bit position values transmitted to, and received from, the controller must be bit/byte-ordered in “Little Endian” format. This means that the least significant bit/byte is last (last to send and last to receive). Byte-order reversal may be required on some platforms. Microsoft Windows, Intel-based Apple Macintosh systems running Mac OS X, and some Intel/AMD processor based Linux distributions handle byte storage in Little-Endian byte order so byte reordering is not necessary before converting to/from 32-bit “long” values. LabVIEW always handles “byte strings” in “Big Endian” byte order irrespective of operating system and CPU, requiring that the four bytes containing a microsteps value be reverse-ordered before/after conversion to/from a multibyte type value (I32, U32, etc.). MATLAB automatically adjusts the endianness of multibyte storage entities to that of the system on which it is running, so explicit byte reordering is generally unnecessary unless the underlying platform is Big Endian. If your development platform does not have built-in Little/Big Endian conversion functions, bit reordering can be accomplished by first swapping positions of the two bytes in each 16-bit half of the 32-bit value, and then swap positions of the two halves. This method efficiently and quickly changes the bit ordering of any multibyte value between the two Endian formats (if Big Endian, it becomes Little Endian, and if Little Endian, it becomes then Big Endian).
6. “Move” commands might have short to long distances of travel. If not polling for return data, an appropriate delay should be inserted between the sending of the command sequence and reception of return data so that the next command is sent only after the move is complete. This delay can be auto-calculated by determining the distance of travel (difference between current and target positions) and rate of travel. This delay is not needed if polling for return data. In either case, however, an appropriate timeout must be set for the reception of data so that the I/O does not time out before the move is made and/or the delay expires.
7. Actual speed for the “Straight-Line Move ‘S’” command can be determined with the following formula: $(5000 / 16) * (\text{sp} + 1)$, where 5000 is the maximum speed in microns/second and “sp” is the speed level 0 (slowest) through 15 (fastest). For mm/second or microns/millisecond, multiply result by 0.001.

NOTES:

INDEX

<i>I</i>	
[SPEED]/LOCK button	22
1	
12.5mm Y Axis	17
2	
25mm X Axis	17
25mm Y Axis	17
5	
50mm X Axis	17
A	
accessories	33
C	
cleaning.....	31
Configuration	14
Configuration switches	15
controller	
cable specs	35
Controls	
TRIO MP-245	14
Power switch.....	14
D	
disclaimer	3
E	
Electrical Connections	13
External control	
Move to specified.....	27
External control	
Move in straight line to specified position at specified speed ('S') command.....	28
External control	
Interrupt straight-line move ('^C') command	28
External control	
Move to specified X-axis position command...	28
External control	
Move to specified Y-axis position command...	29
External control	
Move to specified Z-axis position command...	29
External control	
Setting the angle ('A') command.....	29
External control	
notes.....	29
External-control.....	23
axis position command parameters	25
command sequence formatting.....	24
microsteps and microns (micrometers)	25
protocol and handshaking	24
travel distances and bounds	25
virtual COM port (VCP) serial port settings..	24
External-control commands	
Move to specified.....	27
F	
fuses, mains.....	34
fuses, replacement	
mains.....	3
G	
glassware	
precautions.....	4
H	
headstage Exchange	23
HOME button (while moving to Home).....	22
I	
input	
voltage.....	34
Installation	
electrical connections	13
General	11
headstage mounting.....	13
initial operating instructions.....	13
other accessories	13
ROE/Controller rear panel controls and configuration	14
configuration Switch 9	
electromechanical type.....	17
configuration switches	15
locking out Y-movement homing	15
no calibration homing on power on (FW v2.4+)	16
PULSE button definition.....	16
Pulse mode	16
sensor test.....	16
sensor test (FW < v2.2).....	15
setting axis directionality	15
Switches 7 - 10	16
of X-axis travel length (firmware v2.4+)....	17
power switch	14
Y-axis travel length (FW v2.4+).....	17
TRIO MP-245/M mounting.....	11

Installation	11	Q	
Introduction		quick reference.....	36
Components	10	R	
Overview		Reconfiguration	31
Description	10	changing rotary knob functions on the ROE. 31	
Features.....	10	S	
Overview.....	10	safety warnings	
Introduction	9	mains fuse.....	3
L		safety warnings & precautions	
Locking out Y-movement homing.....	15	operational	3
M		SAFETY WARNINGS & PRECAUTIONS	3
mains		electrical	3
fuses	3, 34	Sensor test.....	16
voltage.....	34	Sensor test (FW < v2.2).....	15
Maintenance.....	31	Setting axis directionality	15
Micropipette Exchange	23	SPEED button	22
Mounting		Switches 7 - 10	16
headstage.....	13	T	
N		technical specifications.....	34
No calibration homing on power on	16	electrical.....	34
notes		V	
user.....	42	voltage	
O		input.....	34
Operations		mains.....	34
control operations	19	W	
headstage exchange	23	warranty	32
maximum positive position values	19	WORK button (while moving to Work).....	22
micropipette exchange	23	X	
mode indications	22	X-axis travel length	
movement knobs disabling and LOCK Mode		25mm	17
.....	22	50mm	17
moving to the Home Position	20	X-axis travel length (firmware v2.4+).....	17
moving to the Work Position.....	21	Y	
pausing Home movements	22	Y-axis travel length	
pausing Work movements	22	12.5mm	17
Pulse Mode and diagonal movement	23	Y-Axis travel length	
ROE axis knob movement speed control....	22	25mm	17
setting Absolute/Relative coordinates mode		Y-axis travel length (FW v2.4+)	17
.....	21		
setting position for HOME and WORK.....	20		
setting the angle of the pipette/headstage			
holder	20		
display	19		
initial startup	19		