# TOPCOM®

## Ear Thermometer 200

**USER GUIDE / HANDLEIDING / MANUEL D'UTILISATEUR /  
BEDIENUNGSANLEITUNG / MANUAL DE USUARIO /  
BRUKSANVISNING / BRUGERVEJLEDNING /  
ANVÄNDARHANDBOK / KÄYTTÖOHJE /  
MANUALE D'USO / MANUAL DO UTILIZADOR / РУКОВОДСТВО  
НА ПОТРЕБИТЕЛЯ / UŽIVATELSKÁ PŘÍRUČKA / ΟΔΗΓΙΕΣ  
ΧΡΗΣΗΣ / HASZNÁLATI ÚTMUTATÓ / INSTRUKCJA OBSŁUGI /  
MANUAL DE INSTALARE SI UTILIZARE / РУКОВОДСТВО  
ПОЛЬЗОВАТЕЛЯ / UŽÍVATEĽSKÝ MANUÁL / KULLANICI  
KILAVUZU**

V 1.0

**CE** 0123

**Manufacturer:** Shenzhen Dongdixin Technology Co., Ltd  
8F-11, No. 35, Hsin-TAI Road  
Chubai City, Hsinchu 302  
TAIWAN

**Distributor:** Topcom Europe NV,  
Grauwmeer 17  
3001 Heverlee, Belgium

This product is in compliance with the essential requirements and other relevant provisions of the directive 93/42/EEC. The declaration of conformity can be found on:  
<http://www.topcom.net/support/cedeclarations.php>

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## 1 INTRODUCTION

Congratulations on the purchase of the Ear Thermometer 200. We hope you will enjoy using the thermometer and the benefits it offers.

The Ear Thermometer 200 uses infrared technology to measure infrared energy emitted from the eardrum and surrounding tissue, and then converts it into temperature value.


## 2 ABOUT BODY TEMPERATURE

Normal human body temperature fluctuates throughout the day; generally, lowest in the morning, highest in the evening.

Different parts of human body also have different temperatures. Therefore, the standards of measuring oral, axillary, rectal and ear temperatures are different. The human body temperature is controlled and adjusted by the hypothalamus in the brain. As the veins flow through the hypothalamus, they also flow through the ear-drum in the mid ear. Thus, measuring the ear temperature (the temperature of the ear drum membrane) can best reflect the body temperature.

Normal body temperature is a range, not a fixed value. When using an IR ear thermometer to take temperature, the normal body temperature range is 35.8 °C - 38.0 °C. Also, people's normal temperature ranges tend to decrease with their age.

However, the normal range is also different from person to person. To determine your and your family members' normal temperature ranges is very important. It can be done easily by using our Infrared Ear Thermometer. Please practice to take temperatures on yourself and your healthy family members to determine their normal temperature ranges.
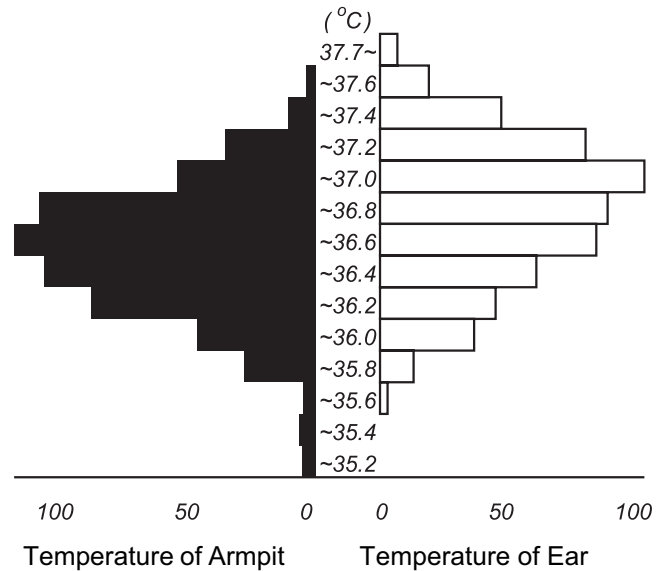
## 3 SAFETY NOTES

- Please read carefully through the following information concerning safety and proper use. Make yourself familiar with all the functions of the equipment. Be careful to keep these advice notes and if necessary pass them on to a third party.
- This thermometer is intended for household use only. The use of this thermometer is not intended as a substitute for physician's consultation.
- When the temperature is higher than 38°, means that the patient may have fever. Contact your doctor.
- Do use a new, clean probe cover while measuring, otherwise the result of measurement will be inaccurate.
- Do not use the thermometer outside the operating ambient temperature of 10°C ~40°C or ambient humidity up over 95% RH.
- For proper measurement, condition the thermometer for 30 minutes in measurement site till room temperature.
- Do not use the thermometer when the ear canal is wet, or when the person is suffering an ear disease.
- Earwax or obstructions will affect the accuracy of the thermometer.
- Sideway sleeping gesture will cause that the ear temperature will rise.
- The thermometer is calibrated and does not need re-calibration if used according to the instructions.
- The measured temperatures may be different in the different ear side. Measure the same ear to ensure an accurate measurement.
- After a shower or exercise, the body temperature will rise. Please take a rest before taking temperature.
- Do not insert the probe into the ear canal with excessive force.
- Don't touch or blow the infrared sensor
- Keep the unit out of children's reach. A self-measurement forced by children may cause ear damage.

## 4 CARE AND MAINTENANCE


- Do not shock, vibrate or drop the thermometer.
- Do not allow any liquid in the thermometer.
- Do not disassemble, repair or remodel the thermometer.
- Keep out of reach of children.
- Clean the main body of the thermometer with a soft and dry cloth.
- Do not use benzene, gasoline, thinner or any other cleaner containing abrasive agents.
- The tip of the probe is the most delicate part of the thermometer, accurate readings from the thermometer can be obtained only if the tip of the probe is clean and intact.
- To clean the tip of the probe, hold the thermometer downwards so that no liquid can enter. Use a cotton swab moistened with alcohol to wipe gently and carefully on the surface of the tip of the probe. Allow an hour drying time before operating the thermometer again.

## 5 BATTERY INSTALLATION

1. Take off the screw and remove the battery cover.
2. Install 2 AAA (1,5V) non-rechargeable batteries following the polarity as shown.
3. Reposition the battery cover and fix the screw


Replace the battery when '  ' is displayed.  
Take out the battery if you do not use it for a long time.  
You have to dispose of the batteries in an environment friendly manner according to your country regulations.
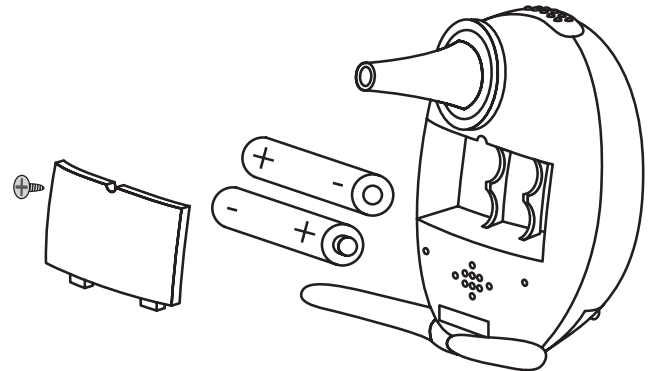
## 6 BUTTONS

1. Start - button
2. Display
3. Talk - button
4. Memory - button
5. °C / °F switch
6. Speaker
7. Reset - button
8. Battery compartment
9. Probe + cover
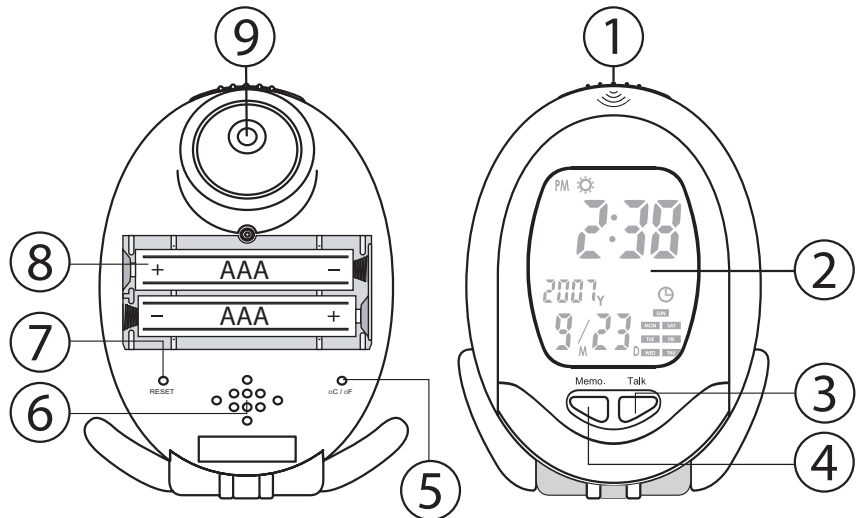
## 7 IMPORTANT


- Replace any low-voltage battery to ensure a full power supply.
- Condition the thermometer for 30 minutes in the measurement site till it reaches room temperature.
- An unexpected fluctuation of ambient temperature will influence the measurement accuracy. When a measurement site shows different temperature with the area where a thermometer is stored, the result will be incorrect.
- Do not perform a temperature measurement in front of an air conditioner.
- Be sure to have a clean probe.
- Be sure to keep the ear clean.
- Measurement shall be done at a rest time to ensure an accurate result. Temperature is likely to increase immediately after exercise or bathing.

## 8 HOW TO USE

- Please make sure to use a new and clean probe cover to ensure accurate measurement.
- Press the **Start** - button (1) to turn the power on. The display will show all symbols.

### Before inserting the probe:

- For children < 1 year : Pull the ear down and back.
  - For children > 1 year to adult : Pull the ear up and back.
- While holding the ear, insert the probe to seal the ear canal, then press the **Start** - button (1) once. You will hear a short beep.
  - When you hear a short melody, the measurement is finished. Take out the thermometer from the ear canal. The measured temperature is displayed. The result will be stored automatically into the memory. If the memory is full (30 records), the oldest measurement result will be removed.

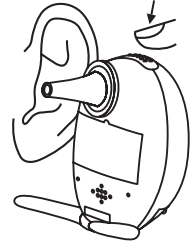It is recommended that you take three temperatures in the same ear under the following circumstances. If they are different, use the highest reading.

- Infants in the first 3 months of life.
- Children under three years old who have a condition such as a compromised immune system and for whom the presence or absence of fever is critical.
- When you are first learning to use the ear thermometer until you are comfortable with the technique and are gaining consistent readings

Important: As with any type of thermometer, slight temperature variations (+/- 0.2 - 0.3°C) can happen, if consecutive measurements are taken.

## 9 DISPLAY ERROR MESSAGES

- **'HI'** The measured temperature is higher than 42,9 °C.
- **'LO'** The measured temperature lower than 32 °C
- **'ErrE'** The thermometer temperature is not within the operation specifications (<10°C or >40°C)
- **'ErrU'** Don't measure after pressing 'Measure'. Follow the steps mentioned in the user guide.
- **'ErrH'** An error is detected during self-check. Remove the batteries for 10 seconds and try again.
- **'

## 10 TIME SETTINGS**

- Power on the unit by pressing the **Start** - button (1) once. You will hear a short beep. Time mode is displayed.
- Press and hold the **Memo** - button (4). The hour digits are blinking.
- Press the **Memo** - button (4) repeatedly to increase the digit.
- Press the **Talk** - button (3) to confirm. The minutes digits are blinking.
- Repeat previous steps for setting the minutes, year, month and day.
- After setting the day, the unit will enter time mode.

## 11 TALK MODE

The result of each measurement is broadcasted through the speaker.  
To broadcast the time, press the **Talk** - button (3) in Time mode.

- To disable the talk mode, press and hold the **Talk** - button (3) for 2 seconds. You hear one beep tone.
- To enable the talk mode, press and hold the **Talk** - button (3) again for 2 seconds. You hear a double beep tone.

## 12 MEMORY

The Ear Thermometer will automatically store the last 30 measurements.

- Power on the unit by press the **Start** - button (1) once. You will hear a short beep. Time mode is displayed.
- Press the **Memo** - button (4). The last result is displayed.
- Press the **Memo** - button (4) repeatly, the display will show the temperature of the previous measurements.
- Press and hold the **Memo** - button (4) repeatly to scroll rapidly.

### 13 TEMPERATURE ENTITY

- Power on the unit by press the **Start** - button (1) once. You will hear a short beep. Time mode is displayed.
- Use a small tool to push the **°C / °F switch** (5) on the back side of the unit.
- A single beep confirms that you have set the temperature to °F
- A double beep confirms that you have set the temperature to °C

### 14 RESET

- Power on the unit by press the **Start** - button (1) once. You will hear a short beep. Time mode is displayed.
- Use a small tool to push the **reset switch** (7) on the back side of the unit.
- All segment will be displayed and all stored information and settings are reset.

### 15 TECHNICAL SPECIFICATIONS

Model	HV-T26
Sensor type	Infrared
Displayed temperature range	32°C ~ 42,9°C (89,6°F ~ 109,3°F)
Operating ambient temperature range	10°C ~ 40°C (50,0°F ~ 104,0°F)
Operating ambient humidity range	30 ~ 80% RH
Storage environment temperature range	-20°C ~ 40°C (-4,0°F ~ 104,0°F)
Storage environment humidity range	30 ~ 85% RH
Display accuracy	0,1°C or 0,1°F
Measurement accuracy	36°C ~ 39°C (+/- 0.2°C) 32C to 35.9 - +/- 0.3 C 39C to 42.9 - +/- 0.3C 96,7°F ~ 102,2°F (+/- 0,4°F) 89,6°F ~ 96,6°F (+/- 0,5°F) 102,3°F ~ 109,3°F (+/- 0,5°F)
Memory	Last 30 measurements
Battery Type	2 x AAA Non- rechargeable - DC 3V
Automatic switch-off	1 Minute after last measurement
Power consumption	< 0,05 (for voice prompt)
Weight	102g (with battery)
Classification 	<ul style="list-style-type: none"> <li>• Class IIa</li> <li>• Internally powered equipment</li> <li>• IPX0</li> <li>• Not suitable for use in presence of flammable anaesthetic mixture with air or with oxygen or nitrous oxide</li> <li>• Continuous operation with short-time loading</li> </ul>

### 16 TOPCOM WARRANTY

#### 16.1 WARRANTY PERIOD

The Topcom units have a 24-month warranty period. The warranty period starts on the day the new unit is purchased. Consumables or defects causing a negligible effect on operation or value of the equipment are not covered. The warranty has to be proven by presentation of the original purchase receipt, on which the date of purchase and the unit-model are indicated.

#### 16.2 WARRANTY HANDLING

A faulty unit needs to be returned to a Topcom service centre including a valid purchase note. If the unit develops a fault during the warranty period, Topcom or its officially appointed service centre will repair any defects caused by material or manufacturing faults free of charge. Topcom will at its discretion fulfil its warranty obligations by either repairing or exchanging the faulty units or parts of the faulty units. In case of replacement, colour and model can be different from the original purchased unit. The initial purchase date shall determine the start of the warranty period. The warranty period is not extended if the unit is exchanged or repaired by Topcom or its appointed service centres.

#### 16.3 WARRANTY EXCLUSIONS

Damage or defects caused by incorrect treatment or operation and damage resulting from use of non-original parts or accessories not recommended by Topcom are not covered by the warranty. The warranty does not cover damage caused by outside factors, such as lightning, water and fire, nor any damage caused during transportation. No warranty can be claimed if the serial number on the units has been changed, removed or rendered illegible.

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