

XI. User registration Cloud Version


1. "INSMART Health" support cloud data storage, if you don't want to lost your data,

register to log in to "INSMART Health" and agree to disclose all information about

IX. Data Comparison

Slide down and click Data Comparison on the home page, select more then two records (Max. 9 records) to compare of the full body composition of their changes.

Click icon on the right top corner to save the comparison records in phone album.

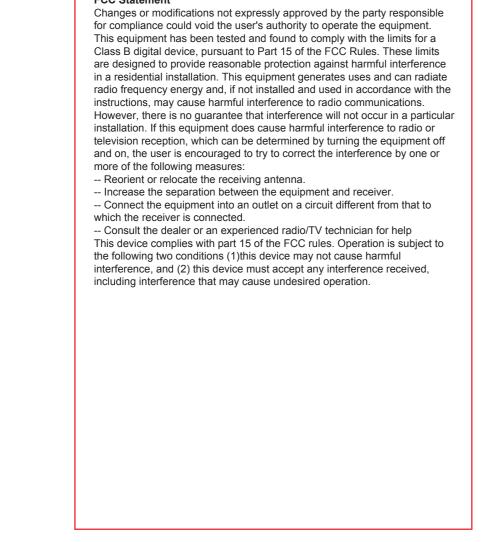

APP FAQ
FREQUENTLY ASKED QUESTIONS
How can I recalibrate the scale?
-Simply step on the scale to turn the power on, step off and let the reading go back to 0.0kg/0.0lb, the calibration will be done automatically. Why doesn't the scale work? Why doesn't the display show data?
Why does the data on the screen disappear in flash?
-Please check whether the insulation piece is taken out, batteries are installed or batteries are in low and need replacing. Display indicates "Err" when measuring body fat:
-Your feet or scale are too wet.
-Overload the maximum capacity(396lb/180kg)
-Body parameter settings should be with Height (3'3 - 7\*37 / 100 -220cm), Aged (10 -99 years old). No body fat data measured when weighing.

-Make sure there are no other Bluetooth connections are turned on (other than the device you are using) within proximity.

-Body fat scale requires weighing with dry and bare feet.

-Please make sure you are stepping on the electro piece or ITO conductive glass of the scale with bare feet the scale with bare feet. -Ensure the Bluetooth function is open. -Please input your personal information with Height(3'3"-7'37 / 100-220cm), Aged (10-99 years old). The scale cannot connect with APP.

-Please check your smart phone model. The App supports Android 4.3 or above with Bluetooth

4.0 for Android devices and IOS8.0 or above for IOS devices.

-For Android 5 or above, please turn on GPS location on the phone before Bluetooth Bluetooth
Connection (the Access Location Info in the APP shows "Allow").
- For IOS11, the Bluetooth is on in setting interface, but actually it may be off in controlling interface. Please make sure the Bluetooth turn on for both setting interface and controlling interface.
- If your phone Bluetooth is connecting to other devices, you need to disconnect them before connect to the scale.

ユーザーズガイド モデル: 260H

ご利用する前に、本説明書をちゃんとお読みください。