
Nokia Body

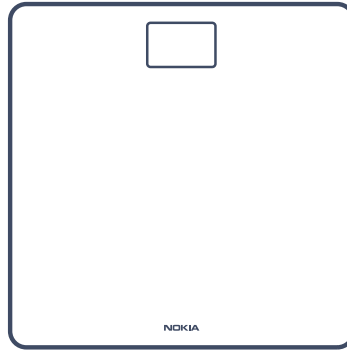QUICK INSTALLATION GUIDE

GUIDE D'INSTALLATION | INSTALLATIONSANLEITUNG | GUÍA DE INSTALACIÓN | GUIDA ALL'INSTALLAZIONE

EN | FR | DE | ES | IT


THANK YOU FOR CHOOSING THE NOKIA BODY

MERCI D'AVOIR CHOISI NOKIA BODY
WIR BEDANKEN UNS FÜR IHRE WAHL ZUM NOKIA BODY
TE AGRADECEMOS QUE HAYAS ELEGIDO EL NOKIA BODY
GRAZIE PER AVER SCELTO LO NOKIA BODY

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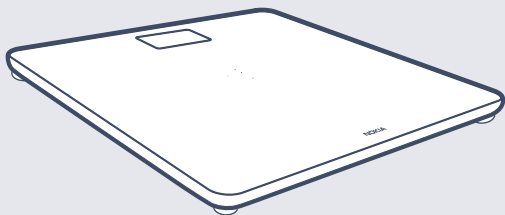
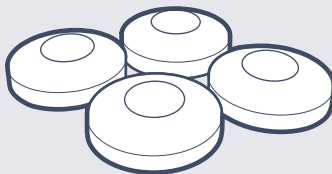
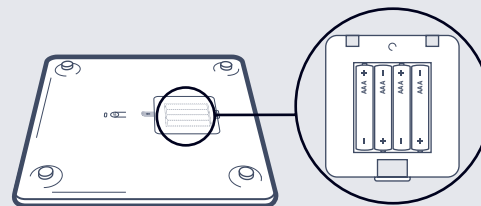
Installation of carpet feet (optional)

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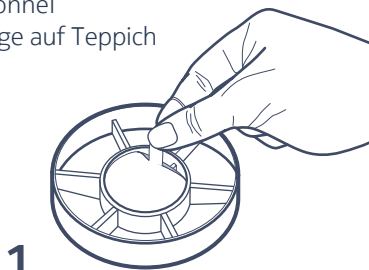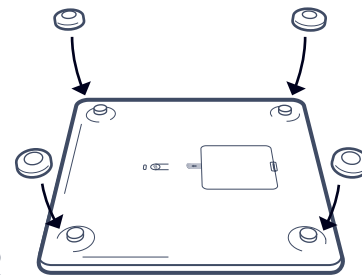
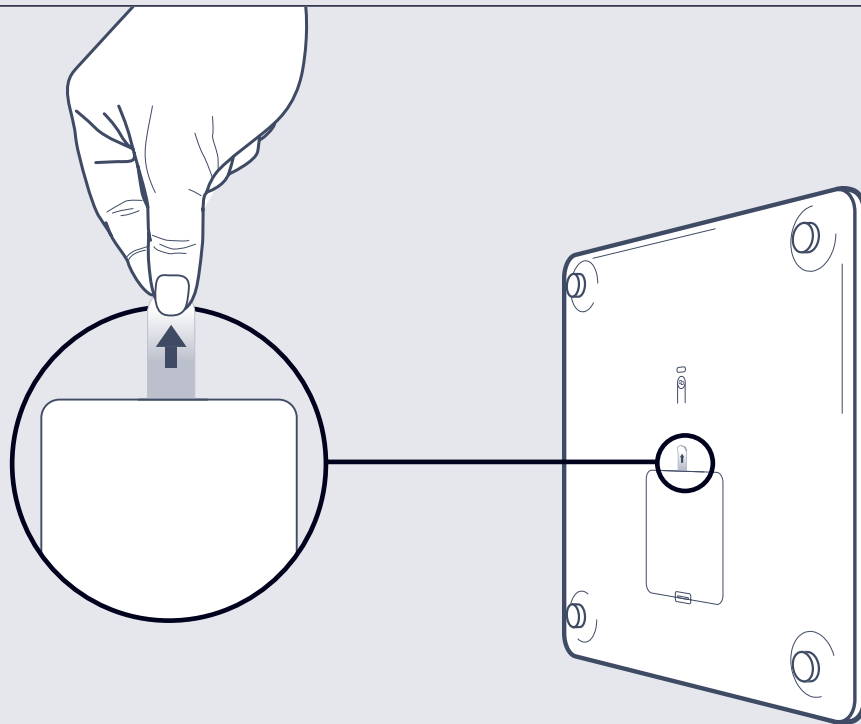
2


Visit healthapp.nokia.com to get the app and follow the instructions.

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
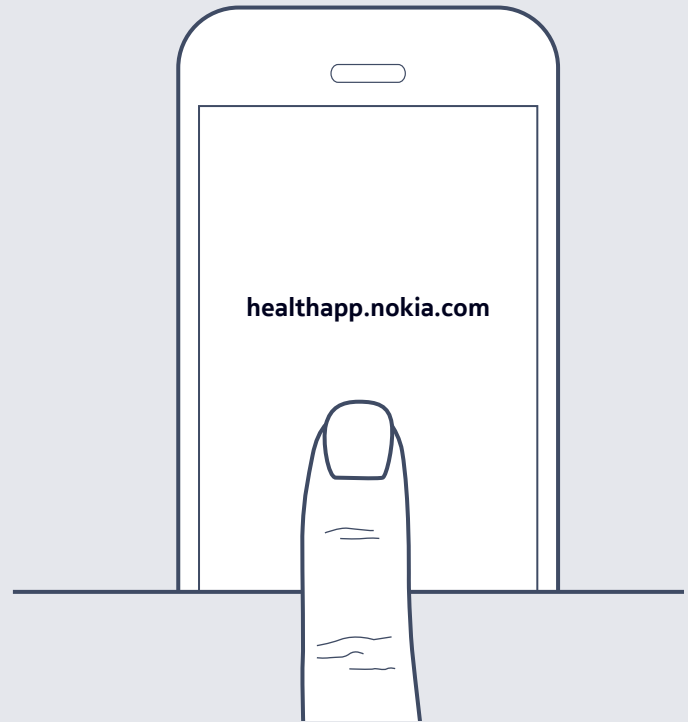
Analyze your Body Mass Index and be notified of Rain and Weather Forecast for the day, while weighing.

- FR** Analysez votre IMC et soyez averti en cas d'alerte pluie ou météo, pendant la pesée.
- DE** Analysieren Sie Ihre BMI und erhalten Sie den Regenbericht für heute beim Wiegen.
- ES** Analiza tu IMC y recibe el Pronóstico de Lluvia para ese día, todo ello mientras te pesas.
- IT** Analizza il tuo BMI e la qualità dell'aria e riceverai, mentre ti pesi, la notifica delle Previsioni meteo del giorno.

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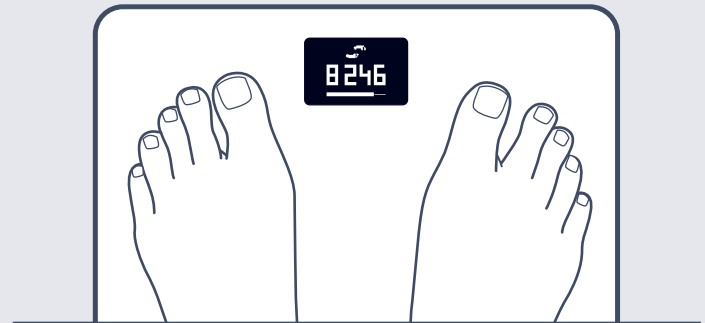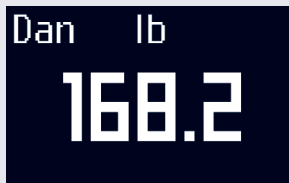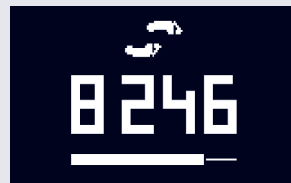
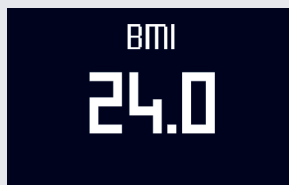
Change the order of the screens and choose what matters most to you.

If you use a Nokia watch, activity tracker or track your steps with the Health Mate™ app, activate the Steps screen to know how much you walked the previous day!

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11

IC: 11411A-WBS06

This device complies with Industry Canada's licence exempt RSSs. Operation is subject to the following two conditions:

- (1) This device may not cause interference; and
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

Radiation Exposure Statement

This equipment complies with Canada radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 20cm between the radiator & your body.

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Garantie | Garantie | Garantia | Garanzia

12

US Nokia™ One (1) Year Limited Warranty - Nokia Body

Withings (Part of Nokia), 2 rue Maurice Hartmann, 92130 Issy-les-Moulineaux, ("Nokia") warrants the Nokia-branded hardware product ("Nokia Product") against defects in materials and workmanship when used normally in accordance with Nokia' published guidelines for a period of ONE (1) YEAR from the date of original retail purchase by the end-user purchaser («Warranty Period»). Nokia' published guidelines include but are not limited to information contained in technical specifications, safety instructions or quick start guide. Nokia does not warrant that the operation of the Nokia Product will be uninterrupted or error-free. Nokia is not responsible for damage arising from failure to follow instructions relating to the Nokia Product's use.

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Notes

Federal Communication Commission Interference Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC Caution: Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

IMPORTANT NOTE:

Radiation Exposure Statement:

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 20cm between the radiator & your body.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Need help?

Besoin d'aide ? | Brauchen Sie Hilfe? | ¿Necesitas ayuda? | Ti serve aiuto?

➔ nokia.com/health/support/body

Mode	Frequency Band (MHz)	Maximum Output Power (dBm)
WLAN	2412-2472	9.87
BT EDR	2400-2483.5	0.85
BT LE	2400-2483.5	9.87

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

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NOKIA

Nokia Body | BMI Wi-Fi Scale

© 2017 Nokia Technologies Ltd. All rights reserved.