





Withings

Made for
iPod iPhone

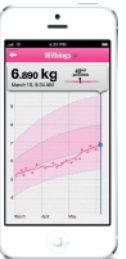

EN / Quick Start Guide for usage with an iOS™ device
FR / Guide d'installation pour utilisation avec terminal iOS™
DE / Installationsanleitung zur Nutzung mit iOS™ Gerät
IT / Istruzioni d'installazione per utilizzo con terminale iOS™
SP / Guía de instalación para utilización con dispositivo iOS™

1

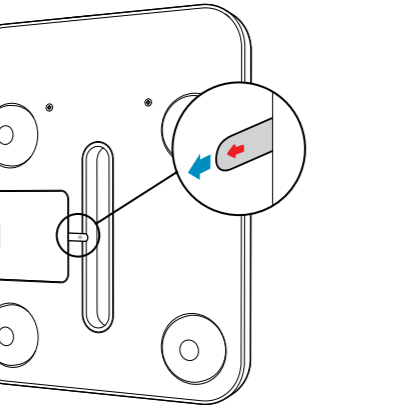


EN / Pull the tab underneath the scale
FR / Retirez la languette à l'arrière de la balance
DE / Ziehen Sie die Lasche unter der Waage ab
IT / Rimuovere la linguetta sotto la bilancia
SP / Retire la pestaña debajo de la báscula

2




EN / Tap the "Settings" app on your iPhone home screen
FR / Touchez l'application "Réglages" dans l'écran d'accueil de votre iPhone
DE / Auf der Startseite Ihres iPhones die App "Einstellungen" berühren
IT / Tocca l'applicazione "Regolazioni" nella schermata iniziale del tuo iPhone
SP / Pulsa la aplicación "Ajustes" en la pantalla de inicio de tu iPhone

3

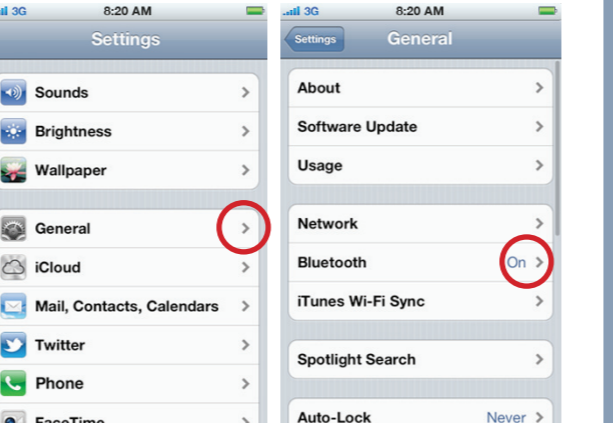


EN / Tap "General", then "Bluetooth"
FR / Touchez "Général", puis "Bluetooth"
DE / "Allgemeine Einstellungen", dann "Bluetooth" auswählen
IT / Tocca "Generale", quindi Bluetooth"
SP / Pulsa "General", y seguidamente "Bluetooth"

4

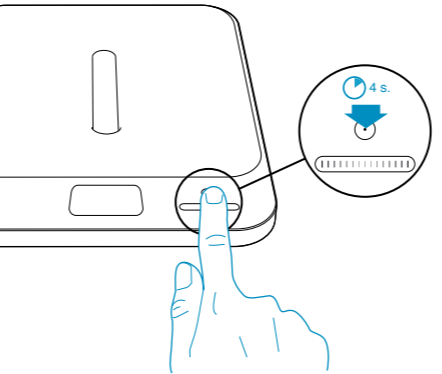

EN / Press the button for 4 seconds
FR / Appuyez 4 secondes sur le bouton
DE / Drücken Sie die Taste für 4 Sekunden
IT / Premere il pulsante per 4 secondi
SP / Pulsa durante 4 segundos el botón

5

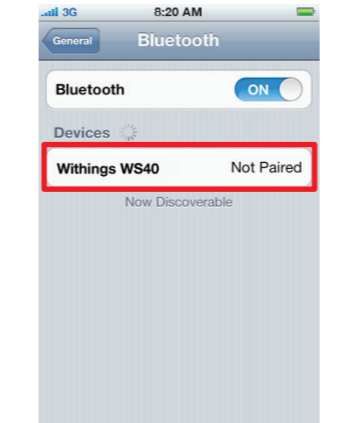


EN / Make sure Bluetooth is turned on, then select the WS40 accessory
FR / Make sure Bluetooth is turned on, then select the WS40 accessory
DE / Sicherstellen, dass Bluetooth aktiviert ist, dann das WS40 auswählen
IT / Comprueba que Bluetooth está activado, y pulsa en el accesorio WS40
SP / Comprueba que Bluetooth está activado, y pulsa en el accesorio WS40

6



EN / Download the Withings app from the App Store
FR / Téléchargez l'application Withings depuis l'App Store
DE / Laden Sie die Withings App aus dem App Store herunter
IT / Scarica l'applicazione Withings da App Store
SP / Descárgate la aplicación Withings desde App Store

7

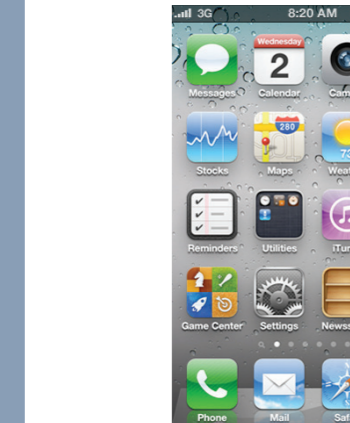

EN / Launch the Withings app and follow instructions
FR / Lancez l'application Withings et suivez les instructions
DE / Starten Sie die Withings App und befolgen Sie die Anleitungen
IT / Lancia l'applicazione Withings et e segui le istruzioni
SP / Inicia la aplicación Withings y sigue las instrucciones

Technical Specifications

Product Specifications

Large graphical screen (2.4" x 1.6")
Touch sensitive interface
Removable baby basket
Four weight sensors
Range: 0 kg to 25 kg (0 - 55 lb)
Graduation: 10 g (0.2 lb)
Display units: kg, lb and oz
Tare function
Dimensions: 27.9 x 23.1 x 2.8 cm (11" x 9.1" x 1.1") without baby basket, 55.9 x 29.7 x 7.4cm (22" x 11.7" x 2.9") with baby basket
Weight: 2.3 kg (5.1 lbs)
Up to 6 profiles with automatic recognition
Wi-Fi connectivity (802.11 b/g/n) compatible with WEP/ WAP/ WPA2-personal security protocols.
Bluetooth® connectivity.


Environmental requirements

Usage conditions:
Temperature: 5°C - 35°C
Relative Humidity: 30% - 80% (without condensation)
Storage conditions:
Temperature: -10°C - +50°C
Relative Humidity: 10% - 90% (without condensation)




EN / Quick Start Guide for usage with an Android™ device
 FR / Guide d'installation pour utilisation avec terminal Android™
 DE / Installationsanleitung zur Nutzung mit Android™ Gerät
 IT / Istruzioni d'installazione per utilizzo con terminale Android™
 SP / Guía de instalación para utilización con dispositivo Android™

1

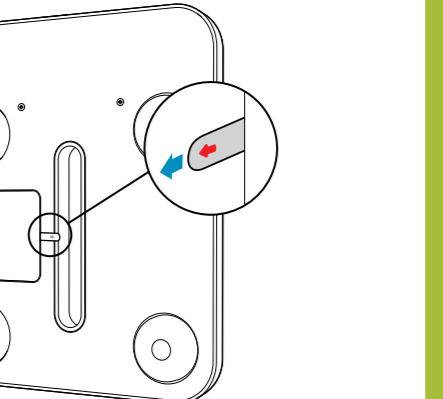

EN / Pull the tab underneath the scale
 FR / Retirez la languette à l'arrière de la balance
 DE / Ziehen Sie die Lasche unter der Waage ab
 IT / Rimuovere la linguetta sotto la bilancia
 SP / Retire la pestaña debajo de la báscula

2

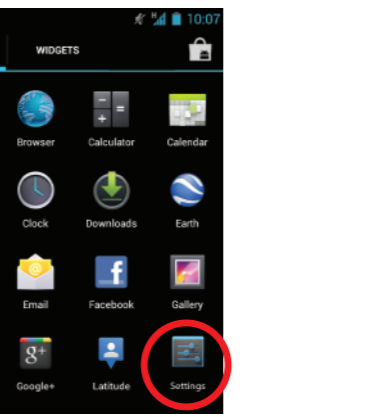


EN / Launch the "Settings" app on your phone
 FR / Lancez l'application "Réglages" de votre téléphone
 DE / Starten Sie die App "Einstellungen" Ihres Telefons
 IT / Lancia l'applicazione "Regolazioni" del tuo telefono
 SP / Inicia la aplicación "Ajustes" de tu teléfono

3




EN / Make sure Bluetooth is turned on and then open the Bluetooth menu
 FR / Vérifiez que le Bluetooth est activé et ouvrez le menu Bluetooth
 DE / Sicherstellen, dass Bluetooth aktiviert ist und dann Bluetooth öffnen
 IT / Verifica che il Bluetooth sia attivato ed aprì il menu Bluetooth
 SP / Comprueba que Bluetooth está activado, y abre el menú Bluetooth

4

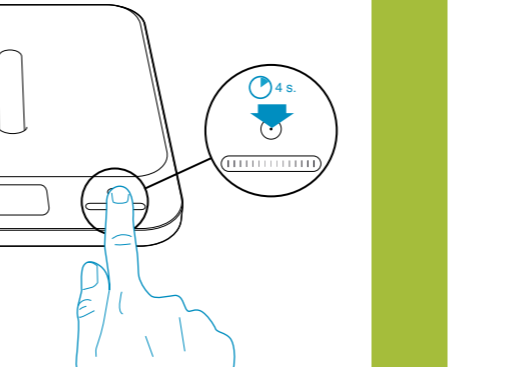


EN / Press the button for 4 seconds
 FR / Appuyez 4 secondes sur le bouton
 DE / Drücken Sie die Taste für 4 Sekunden
 IT / Premere il pulsante per 4 secondi
 SP / Pulsa durante 4 segundos el botón

5

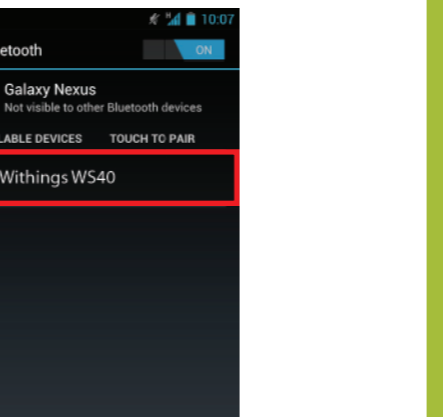


EN / Select the WS40 accessory
 FR / Sélectionnez l'accessoire WS40
 DE / WS40 Zubehör auswählen
 IT / Seleziona l'accessorio WS40
 SP / Selecciona el accesorio WS40

6

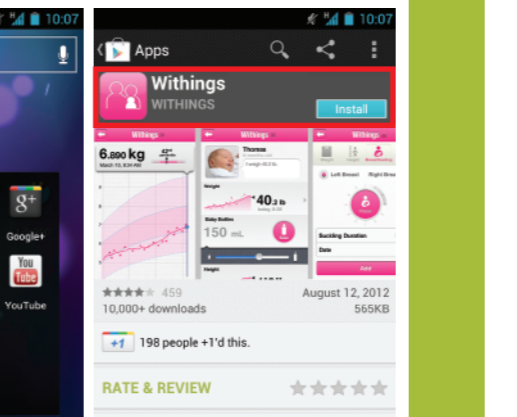

EN / Download the Withings app from Google Play
 FR / Téléchargez l'application Withings sur Google Play
 DE / Laden Sie die Withings App aus Google Play herunter
 IT / Scarica l'applicazione Withings su Google Play
 SP / Descárgate la aplicación Withings en Google Play

7

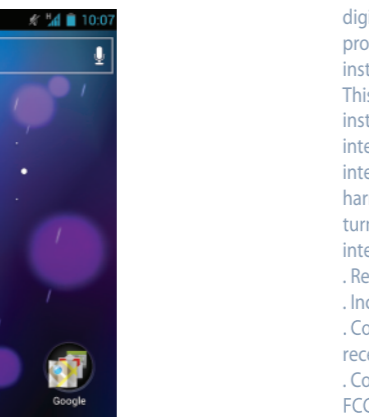

EN / Launch the Withings app and follow instructions
 FR / Lancez l'application Withings et suivez les instructions
 DE / Starten Sie die Withings App und befolgen Sie die Anleitungen
 IT / Lancia l'applicazione Withings et e segui le istruzioni
 SP / Inicia la aplicación Withings y sigue las instrucciones

Federal Communication Commission Interference Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of the following measures:

- . Reorient or relocate the receiving antenna.
- . Increase the separation between the equipment and receiver.
- . Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- . Consult the dealer or an experienced radio/TV technician for help.

FCC Caution: To assure continued compliance, any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment. (Example - use only shielded interface cables when connecting to computer or peripheral devices).

FCC Radiation Exposure Statement

This equipment complies with FCC RF radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with a minimum distance of 20 centimeters between the radiator and your body.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter. The antennas used for this transmitter must be installed to provide a separation distance of at least 20 cm from all persons and must not be co-located or operating in conjunction with any other antenna or transmitter.



FCC ID: XNA-WS-40

Federal Communication Commission Interference Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC Caution: Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

IMPORTANT NOTE:


Radiation Exposure Statement:

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 20cm between the radiator & your body. (note: the radiator is turned off when weight-in measurement occurs)


This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

The availability of some specific channels and/or operational frequency bands are country dependent and are firmware programmed at the factory to match the intended destination.


Baby scale / Pèse-bébé

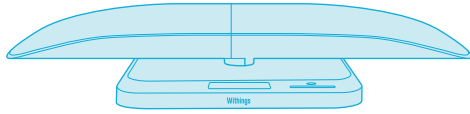


1

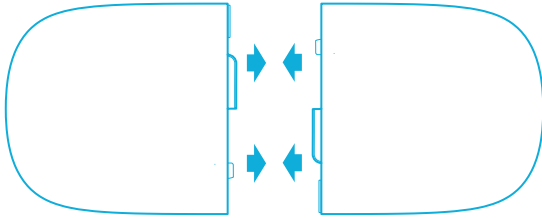

2

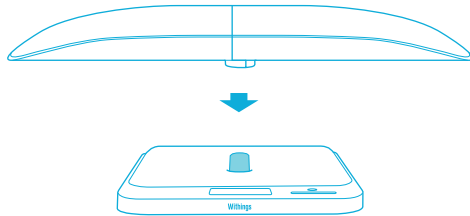


3

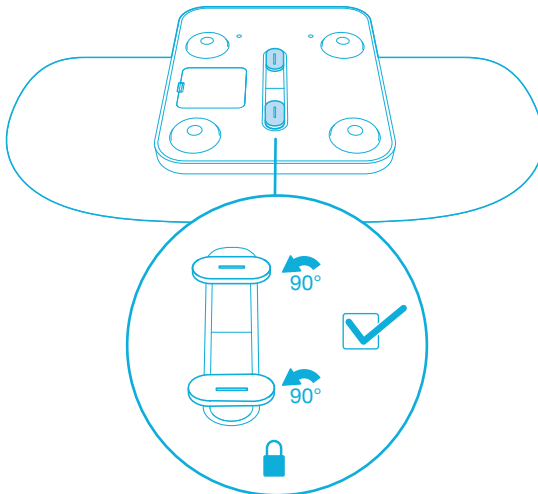

Kid scale / Pèse-enfant



55 lb max / Max 25 kg

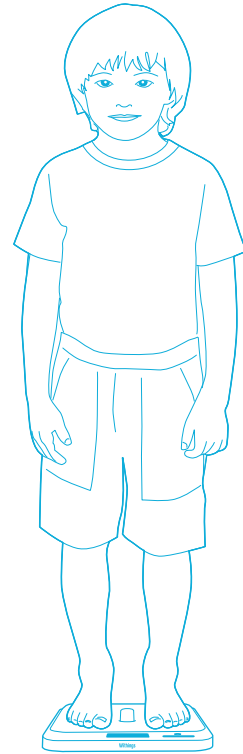

<http://support.withings.com>

Withings