
# Withings

Smart Body Analyzer  
WS-50
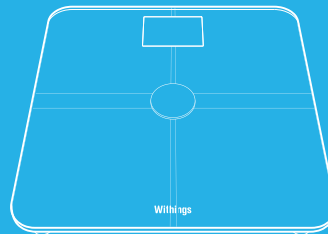
## Quick Installation Guide

---

**FR** Guide d'installation    **ES** Guía de instalación  
**DE** Installationsanleitung    **IT** Guida all'installazione


## NEED HELP

BESOIN D'AIDE | BRAUCHEN SIE HILFE | NECESITAS AYUDA | TI SERVE AIUTO


➔ <http://support.withings.com>

## BOX CONTENTS

CONTENU DE LA BOÎTE | INHALT | CONTENIDO DE LA CAJA | CONTENUTO DELLA CONFEZIONE

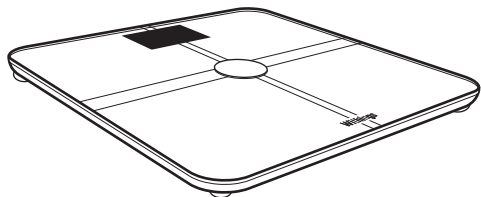
**Smart Body Analyzer  
WS-50**

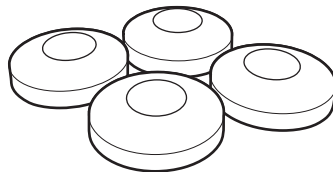**Carpet feet**

FR Pieds moquette ES Pies para moqueta  
DE FüÙe für Teppich IT Piedini per moquette

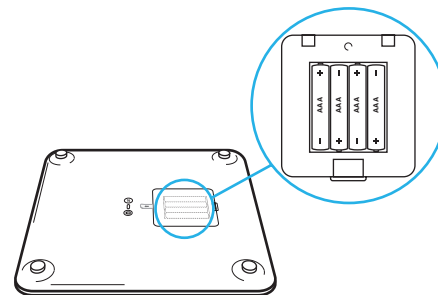**Alkaline cells already installed**

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

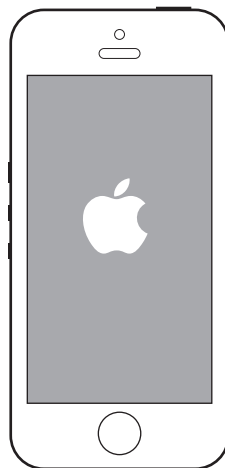
iPhone,  
iPad & iPod touch

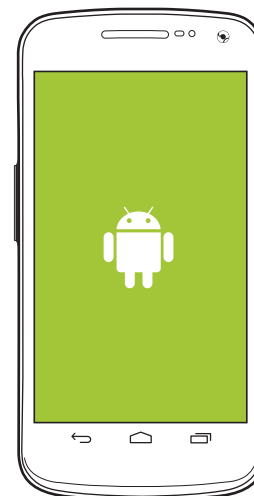Android
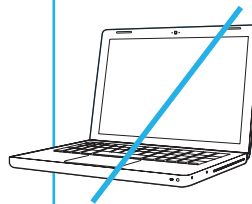
## Installation using your smartphone/tablet only. Computer-free.

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
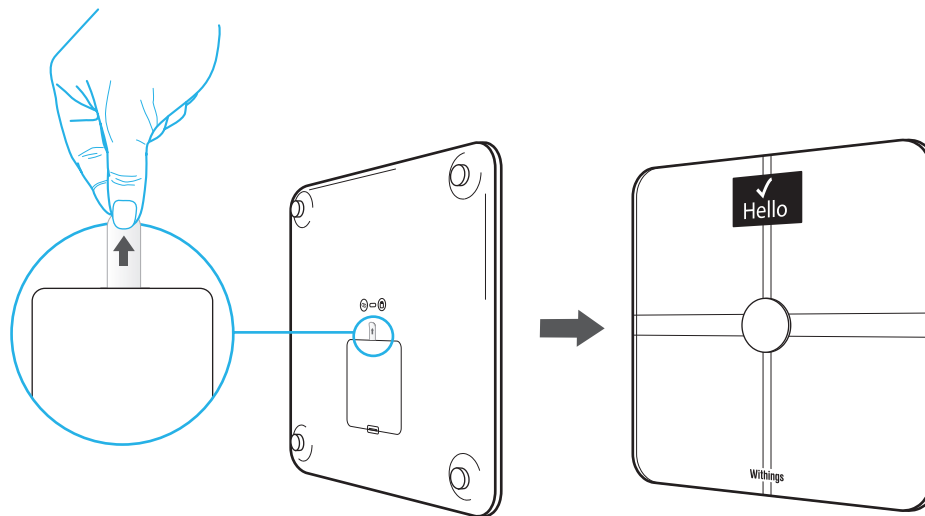
2

## Press the pairing button for two seconds

- FR Appuyez 2 secondes sur le bouton d'appairage  
DE Drücken Sie 2 Sekunden lang die Pairing-Taste  
ES Pulsa durante 2 segundos el botón de aparejar  
IT Premere 2 secondi il tasto di accoppiamento
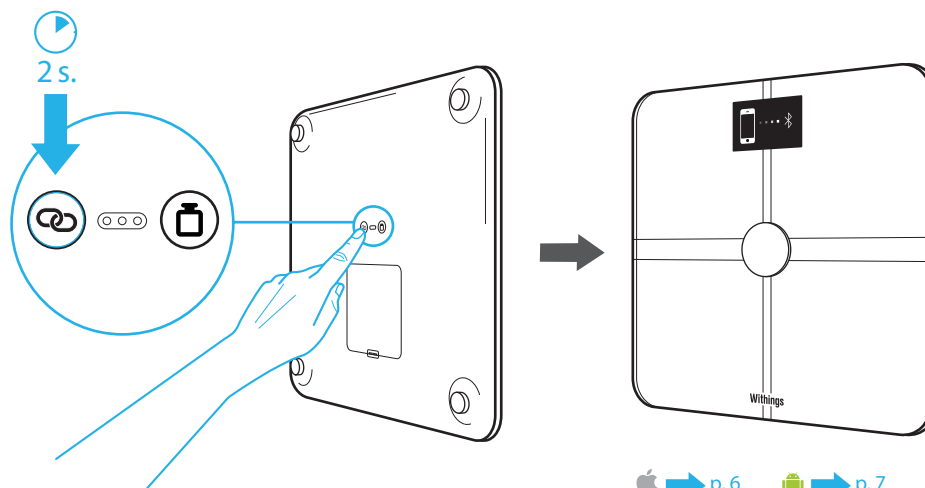
# iPhone iPad, iPod Touch

## 3

### Settings / Bluetooth / ON → "Withings WS50"

- FR Réglages / Bluetooth / ON → "Withings WS50"
- DE Einstellungen / Bluetooth / ON → "Withings WS50"
- ES Ajustes / Bluetooth / ON → "Withings WS50"
- IT Regolazioni / Bluetooth / ON → "Withings WS50"

Settings

## 4

### Install and launch the Withings app from the App Store then follow the instructions

- FR Installez et lancez l'appli Withings, puis suivez les instructions
- DE Installieren und starten Sie die Withings-App und folgen Sie dann den Anweisungen
- ES Instala y abre la aplicación Withings y sigue las instrucciones
- IT Installare e lanciare l'app Withings, poi seguire le istruzioni

# ANDROID


## 3

### Settings / Bluetooth / ON → "Withings WS50"

- FR Réglages / Bluetooth / ON → "Withings WS50"
- DE Einstellungen / Bluetooth / ON → "Withings WS50"
- ES Ajustes / Bluetooth / ON → "Withings WS50"
- IT Regolazioni / Bluetooth / ON → "Withings WS50"



Settings

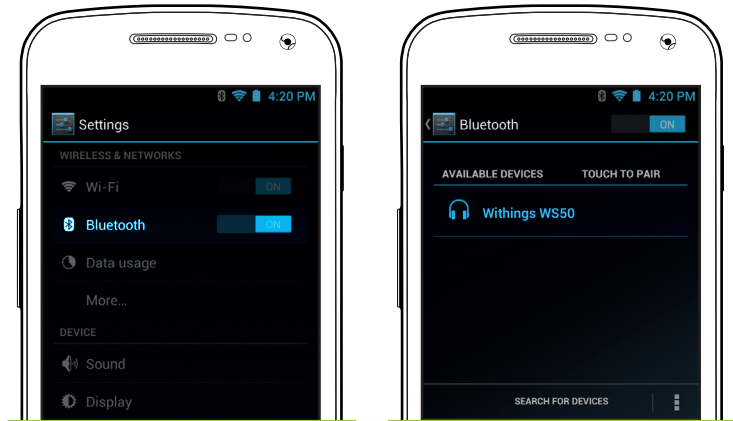
## 4

### Install and launch the Withings app from the Google Play Store then follow the instructions

- FR Installez et lancez l'appli Withings, puis suivez les instructions
- DE Installieren und starten Sie die Withings-App und folgen Sie dann den Anweisungen
- ES Instala y abre la aplicación Withings y sigue las instrucciones
- IT Installare e lanciare l'app Withings, poi seguire le istruzioni

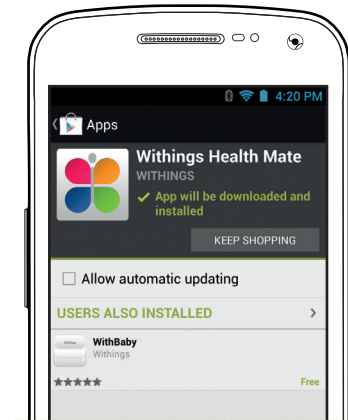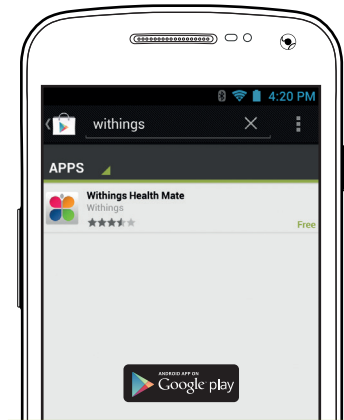Withings
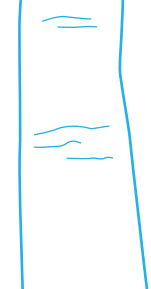
## IMPORTANT ADVICES

CONSEILS IMPORTANTS | WICHTIGE HINWEISE | CONSEJOS IMPORTANTES | SUGGERIMENTI IMPORTANTI

- This product is not suitable for people wearing a pacemaker or other internal medical device.
- The Withings scale data must be considered purely as reference values and not as a medical opinion.

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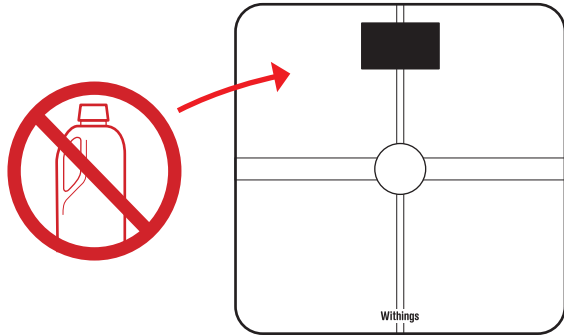
**- Never use detergents to clean your scale. You may irreversibly damage the coating that protects the electrodes used to measure your body analysis.**

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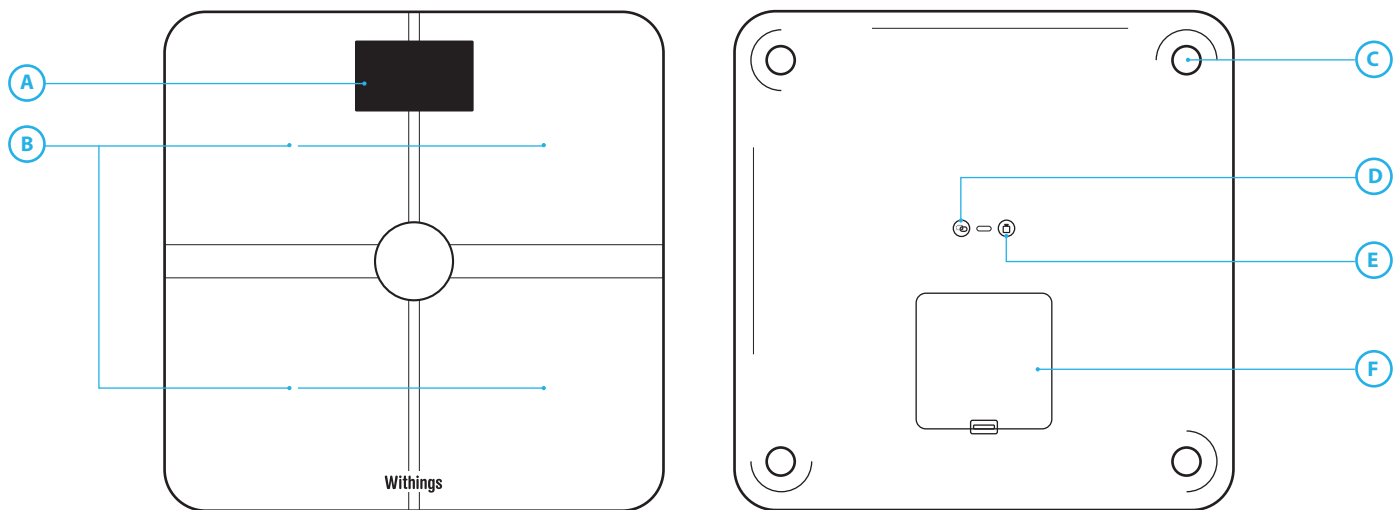
## **FEDERAL COMMUNICATIONS COMMISSION INTERFERENCE STATEMENT**

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

### **CAUTION:**

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

### ***FCC Radiation Exposure Statement***

This equipment complies with FCC RF radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with a minimum distance of 20 centimeters between the radiator and your body.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

The antennas used for this transmitter must be installed to provide a separation distance of at least 20 cm from all persons and must not be co-located or operating in conjunction with any other antenna or transmitter.

# Withings

WS-50 QIG | v1.0 | February 22nd, 2013

©2013 Withings. All right reserved.

iPad, iPhone and iPod touch are trademarks of Apple inc.

Android is a trademark of Google Inc.