YUNMAI mini


User Manual for YUNMAI Smart Scale



Name

For personal health data gathering, you are required to enter your real name when first using the App of YUNMAI Smart Scale.

Click the text box and enter your name.

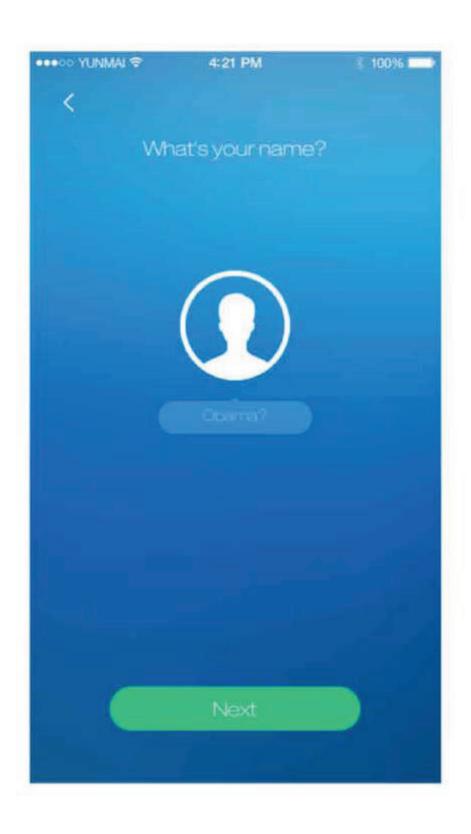

Age

Your age is extremely important to YUNMAI Smart Scale, since most health data is closely related to a person's age. For the sake of your privacy, we will surely keep it private.

Swipe your finger up or down the screen to choose your age.



Height

Your height is extremely important to YUNMAI Smart Scale. Make sure you choose the right one .


Swipe your finger up or down the screen to choose your height.



Physical Examination Results

When you step on the YUNMAI Smart Scale, you can use the App not only to check your body weight, body fat, BMI, etc., but also to receive an overall data of your physical condition.

With a gentle tap on each of the body compositions, you can find their detailed imformation.

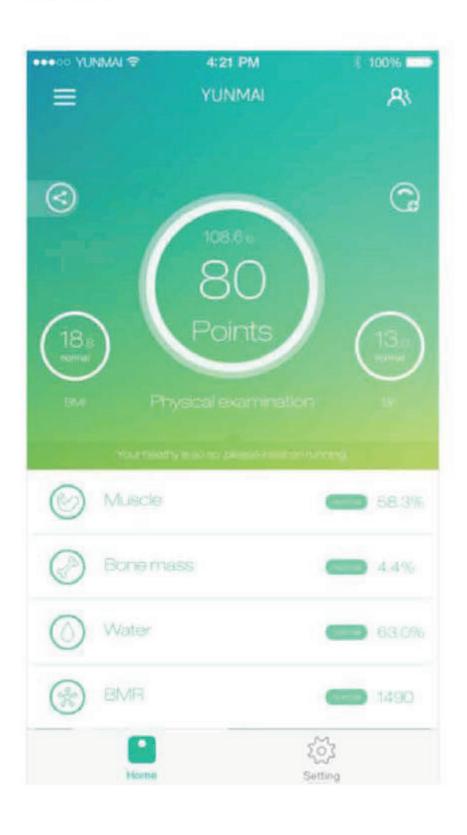

Briefing


Every time after you use YUNMAI Smart Scale, the App will automatically store your date in the cloud. You can easily check body weight or body fat during the past week, month or even the past whole year.

By taping the arrow-shape button, you can see your daily data within a week.





Setting Goals

YUNMAI Smart Scale enables you to set your goal, to be weight losing, weight gaining or weight staying. Enter your target and it's done. Swipe your finger left or right on the screen to set your target body weight.

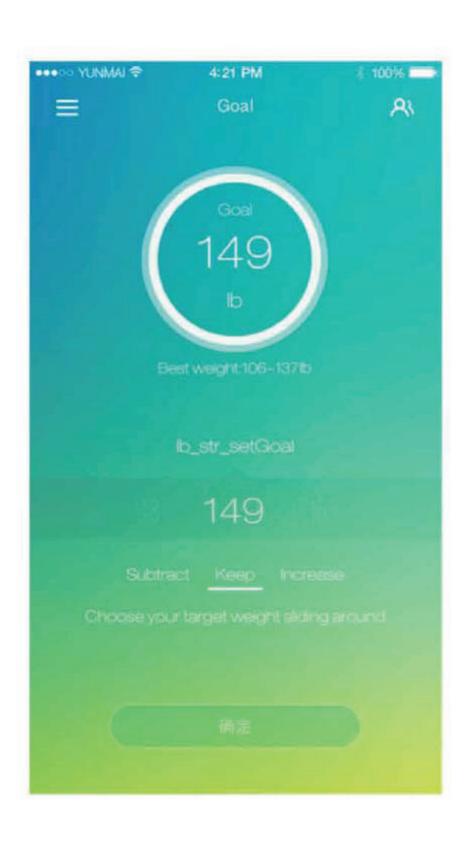

Thank you for purchasing the YUNMAI Smart Scale. Below is the product description.

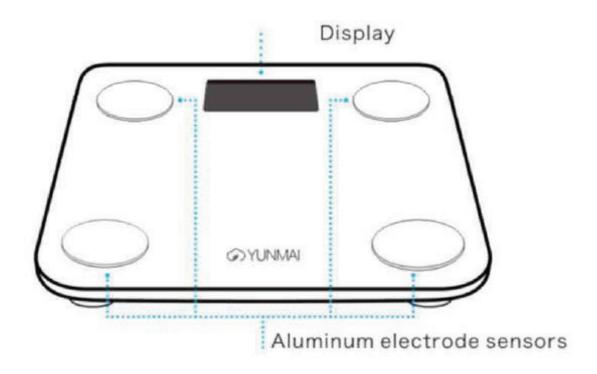


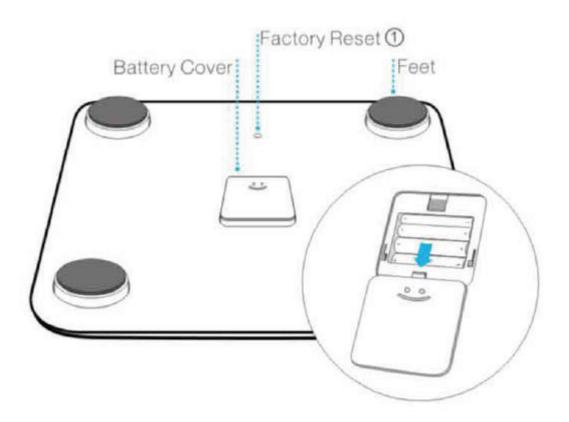

① Important! The factory reset must be used for assembly by staffs only.

Correct Standing Posture

To get accurate results, you must be barefoot and step right on the aluminum electrode sensors.

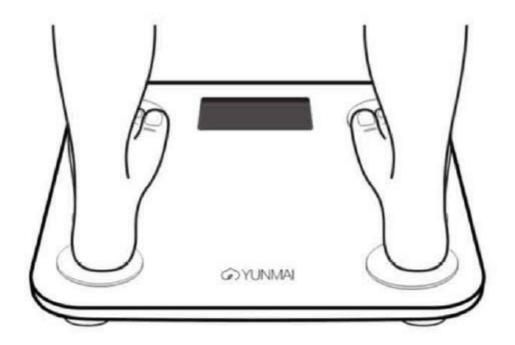

Note: To measure body statistics, you must connect APP.

Product Specifications

Product Name: YUNMAI mini Smart Scale

Model: M1501

Weight Range: 3KG -180KG

Measuring Accuracy: +-50g

Weighing Unit: KG/LB

Working Temperature: 10° C - 40° C

Net Weight: 0.9KG Shell Material: ABS

Power: 4X1.5 V AAA Batteries

Size: 260 mm X 260 mm X 19 mm

Connect with iphone

iPhone4S+, iPhone, iPod Touch, iPad iOS 7.0 and above.

1.Pls click the Setting icon



3. Turn on Bluetooth



2.Choose Bluetooth item



 Click YUNMAI App to start healthy journey



Connect with Android

YUNMAI Smart Scale available for Android 4.3 or above. Pls keep Bluetooth turn on mode to automatically connect with smart scale.

Not support version below Android 4.3.

Important

Use the YUNMAI Smart Scale while barefoot. The scale cannot measure body fat percentage if you're wearing shoes or socks.

YUNMAI Smart Scale should not be used to diagnose or treat any medical situation.

Limitation

The weighing range of YUNMAI Smart Scale: 6.6 lbs (3 kg)-396.8 lbs (180kg) Accuracy: 0.1 kg / 0.2 lb

	Length	Width	Depth	Weight
YUNMAI Smart Scale	260 mm	260 mm	19 mm	0.9 kg
	12.8 in	12.8 in	1.2 in	4.76 lbs

FCC Statement

FCC ID: 2ADEB - M1501

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.
 This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:
- (1) This device may not cause harmful interference, and(2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Read this manual before use. Keep it for future reference

Intended Use: YUNMAI Smart Scale is a digital monitor intended for use in measuring weight and body fat mass in generally healthy adults beyond 18 years old. It is intended for use at home/domestic setting only.

Care and Maintenance

Please clean the surface of the YUNMAI Smart Scale with a damp cloth. Do not use any detergents and cleansers for cleaning the Aluminum alloy electrode plate. Aluminum electrode plate can measure body fat

Read this manual before use. Keep it for future reference

Intended Use: YUNMAI Smart Scale is a digital monitor intended for use in measuring weight and body fat mass in generally healthy adults beyond 18 years old. It is intended for use at home/domestic setting only.

More Information

You can find more information on official website: www.lyunmai.com
Join the community and know more detailed functional information.



Scan the QR code to get this APP Enjoy more functions of YUNMAI

Customer Service

YUNMAI Smart Scale provides permanent technical support and hardware warranty period within one year.

Official: support@iyunmai.com

Thank you for choosing YUNMAI Smart Scale.