


TEO ClearPAGE™ Precast Gel Running Instructions for SDS PAGE and DNA/Native Gels

Table of Contents

General Information	1
Running buffer preparation	1
Sample preparation	2
Gel preparation	2
Loading of samples onto gel	2-3
Running the gel	3
Gel staining	4-5
Gel drying	6-7
Gel Blotting	8
Buffer formulations	9-10
2-Color Marker	11
Ordering Information	12-14

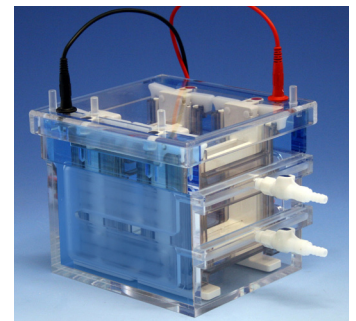


ClearPAGE™ Gel Running Instructions

FOR USE WITH: C.B.S. Scientific's Combo 2-place Dual Cool System (DCX-700) and
Combo 4-place Quadra Systems (QNC-700 & QNX-700)

GENERAL INSTRUCTIONS

- 1) Use **ONLY** ClearPAGE buffers with ClearPAGE TEO gels. Ordering information for these reagents is on pages 12-14. Alternatively, you can make your own by following the buffer formulations on pages 9-10.
- 2) The running and loading buffers must be consistent with the type of gel being used (SDS PAGE or DNA/Native)
- 3) Do not substitute or mix reagents with those from other sources.
- 4) Use only ultrapure water for washing and dilution purposes.
- 5) These products should only be used by qualified personnel who have had laboratory safety instruction.
- 6) The complete instructions should be read and understood before attempting to use the product.
- 7) Products should be used before the "Best by" date on the product label.
- 8) **Wear gloves at all times when handling gels, buffers, etc.**
- 9) Before running first gel in the DCX-700, QNC-700 or QNX-700 please read the Dual Cool/Quadra Instruction manual for important safety and general use instructions. To open DualCool/Quadra, push thumbs down on white/clear posts while pulling lid up with fingers as shown on right. **DO NOT PULL UP ON LEADS!**

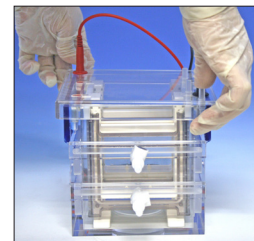

RUNNING BUFFER PREPARATION

Type of Gel	Running Buffer Cat. #'s	Dilution factor
ClearPAGE Neutral pH -TEO	SDS-Non-Reducing: (cat. # FB50053 or FB50500)	20X
	SDS-Reducing: (cat. # FB60053 or FB60500)	20X
	Turbo SDS: (cat.# CB13500)	20X
ClearPAGE "Classics" TGS	Transfer and ClearPAGE "Classics" (FB82500)	20X
DNA/Native	DNA/Native: (cat. # GB61053 or GB61500)	20X

ClearPAGE™ SDS or ClearPAGE DNA/NATIVE Running Buffers **MUST** be used to run these gels. To make 1 liter of 1x running buffer (for example: dilute 50ml 20x Running Buffer with 950ml ultrapure water or dilute 100mls of 10X buffer with 900mls of ultrapure water) for a run. For the ClearPAGE Running Buffer bottle included in the Sampler Kits, dilute by adding its contents to 950ml ultrapure water. Remember, to save time and money you can use reconstituted ClearPAGE running buffer for the bottom reservoir by following instructions below:

RECYCLING BUFFER INSTRUCTIONS: The buffer may be reused on the anode side, but fresh buffer is **ALWAYS** required in the cathode chamber. When reusing buffer make sure used buffer from upper core gets recycled as well. This can be accomplished by removing the gel cassettes while the core is still in the unit and letting the buffer drain into the lower chamber. Pour this reconstituted buffer into a separate container for reuse in the **anode** chamber of your next gel run! This buffer can be used time and time again. For fresh upper core buffer dilute 10 ml of ClearPAGE™ Running Buffer (20x) to 200 ml with ultrapure water.

SAMPLE PREPARATION

Sample Preparation - Prepare samples using ClearPAGE sample buffers for optimal resolution. Other buffers can be used but resolution will suffer. 1 part ClearPAGE 4x LDS sample buffer can be added to 5 parts existing sample to aid with loading and stacking into the gel.

- 1) For DNA/Native Gels, ClearPAGE 4x DNA/Native Sample Buffer (Cat # GB33002) is required for sample preparation.
- 2) For SDS Gels, ClearPAGE 4x LDS Sample Buffer (Cat # FB31002/FB31010) is required for sample preparation.
 - a) Allow Sample Buffer to equilibrate to room temperature
 - b) For **reducing** conditions
 - i) Combine Sample and ultrapure water (if necessary) to achieve 65% of the final volume
 - ii) Add 25% of the final volume of 4x LDS Sample Buffer
 - iii) Add 10% of the final volume of reducing agent (e.g., Cat # FB32001 DTT Reducing Agent 10x).
 - iv) Heat the samples for 10 minutes at 70°C (preferred) or 3 minutes at 100°C.
 - v) Once reduced, run samples within 2 hours to prevent reoxidation.
 - c) For **non-reducing** conditions
 - i) Combine Sample and ultrapure water (if necessary) to achieve 75% of the final volume
 - ii) Add 25% of the final volume of 4x LDS Sample Buffer (Cat # FB31002/FB31010)
 - iii) Heat the samples for 10 minutes at 70°C (preferred) or 3 minutes at 100°C.

Below is an example on how to prepare 100 μ l of 1 mg protein/ml running sample from a 5 mg protein/ml initial sample.

Solution	Reduced	Not Reduced
Sample	20 μ l	20 μ l
Water	45 μ l	55 μ l
Sample Buffer	25 μ l	25 μ l
Reducing Agent	10 μ l	None

Well Type	Max. Volume	40% Max. Vol.
12-well	35 μ l	14 μ l
17-well	17 μ l	6.8 μ l
1 Prep well	400 μ l	280 μ l
2D-well	300 μ l	210 μ l

GEL PREPARATION


- 1) The comb has been removed from these gels prior to shipping and there is also no tape to remove. Just before use remove the gel cassette from its plastic storage bag and shake gel off upside down.
- 2) Rinse entire gel cassettes using ultrapure water with a wash or squeeze bottle.
- 3) Rinse wells two times using ultrapure water with a wash bottle, shake to remove water from wells after each wash.

LOADING SAMPLES ONTO THE GEL:

- 1) To load samples slide precast gel cassette into core assembly so that the notched plate faces towards the core upper reservoir.
- 2) After rinsing wells, fill with ultrapure water. Density differential between water and sample buffer aids protein stacking into the gel. If filling with running buffer, leaving for an extended period will adversely effect resolution.
- 3) When running one gel, use white plastic adaptor plate to seal the other side. Close core doors and press down on white latches. Fill upper core reservoir



with \approx 200mls of buffer. Acrylamide gel well fingers are red (SDS Protein) or green (DNA/Native) and extend above cassette for ease in loading and the wells are numbered. Using thin pipette tips underlay the samples near the bottom of the well as shown at right. Red arrow indicates how acrylamide gel well fingers extend above cassette for ease in loading. Maximum sample capacity is 35 μ l for 12-well gels and 17 μ l for 17-well gels.

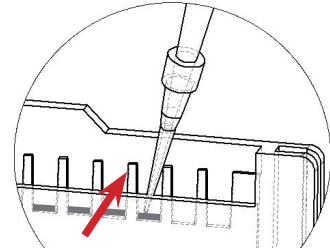

Note: Gel samples may also be loaded on the bench-top and then transferred to the DCX-700. With gels clamped into DCX/Quadra core, fill the wells with running buffer and underlay sample with thin pipette tips prior to transferring to your DCX-700.

RUNNING THE GEL

1) After filling the upper core reservoir with fresh running buffer and loading the samples, add a minimum of 400ml running buffer to the DCX-700 lower reservoir. See chart below if using freezer blocks. To the Quadra Systems add a minimum of 550mls to the lower reservoir.

# of freezer blocks used	DCX-700 lower buffer volume
0	810 mls (max)
1	685 mls (max)
2	560 mls (max)

2) The red or green dye in the loading area of the gel will run ahead of the bromphenol blue. Run the gel(s) until the blue dye front nears the bottom of the cassette as follows:

For DNA/Native Gels, run the gel(s) at 180VDC constant until the blue dye nears the bottom of the cassette (40-80 minutes depending on gel percentage).

For SDS Gels, run the gel(s) until the blue dye front nears the bottom of the cassette as follows:

Run Voltage	Starting Current	Ending Current	Approx. Run Time
180VDC	90mA/gel	40mA/gel	30-75 minutes

AFTER THE RUN




1) Turn off the power supply and disconnect the leads from the power supply.
2) Remove the safety cover from the unit by placing thumbs on the white posts next to the red and black connectors, then pushing down while pulling up with fingers under lid as shown in figure at right. **Do not remove safety cover by pulling up on leads wires!**

3) Leave core in place and pull up on core latches and open gel door. Remove gel cassette(s). Allow upper buffer to mix with lower buffer and **save for reuse as anodal side of next run.**



4) **Completely immerse core and all DCX-700 or Quadra components with deionized or UP water and soak for at least a few minutes.**

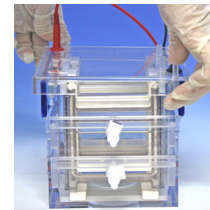

OPENING THE CASSETTE & PREPARING GEL FOR STAINING


1) Using wedge plate separator, dime, comb or gel knife, gently lever apart all four corners of the cassette, first on one side, then the other. Not much force is required to open these cassettes.

GEL STAINING

These composite polymer gels are very strong and change shape minimally when placed in different solutions. However, the composite polymer tends to hold SDS strongly compared to other 1mm thick mini-gels, so a pre-wash or extra staining time is necessary to remove the SDS which competes with the proteins for the stain in more sensitive staining methods. **This step is NOT necessary when using the ClearPAGE™ Instant Blue Stain.**

1) Staining with ClearPAGE™ Instant Blue Stain

- BEFORE USE:** Mix the Instant Blue solution immediately before use by gently inverting the bottle a few times (do not shake the bottle to mix the solution). One may rinse the gel surfaces briefly with ultrapure water, but do NOT wash the gel, as standing in water before staining will reduce band sharpness. Fixing is also NOT recommended.
- GEL STAINING:** Place each gel in a separate small container (lids for yellow pipette tips work well). Use the following amount of stain for different size gels: Use 25ml per gel for gels approximately 8 x 8 cm (from 10 x 10 cm cassettes, such as ClearPAGE or Invitrogen Gels). Use 20 ml per gel for gels approximately 6 x 8cm (from 8 x 10cm cassettes, such as BIO-RAD, Precise or Gradipore gels) Larger Volumes are required for larger gels and for larger containers.
- Shake gently until desired band intensity is achieved. Generally, 50ng or less can be seen in 10 minutes. Band intensity reaches a maximum in 30 to 60 minutes depending on gel thickness.
- Destaining does NOT improve sensitivity, but the gels may be washed in ultrapure water for 15 minutes to remove the free stain from the gel and get clear backgrounds. The gel may be stored in ultrapure water.

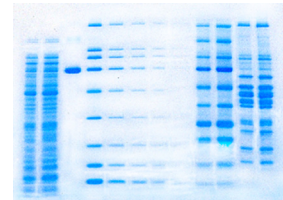

Instant Blue stained ClearPAGE 12% SDS 12-well gel (30 minutes), followed by brief water wash (5 minutes).

Lane	Samples loaded
1-2	E. Coli Extract
3	Bovine Serum Albumin (2µg)
4-8	Broad Range MW Marker 500/200/100/50/25ng per band
9-10	ClearPAGE Two-Color SDS™ Marker 2µl / 4µl
11-12	Chicken Muscle Extract

2) Staining with Coomassie™ Brilliant Blue R-250

When using R-250 stain solutions (40% ethanol/10% acetic acid/0.1% Brilliant Blue R-250), follow standard protocols or use the following recommended by ClearPAGE. Use 20% ethanol/5% acetic acid as the R250 destain solution, half the usual destain concentration, combined with an absorbent paper such as paper towel or lab wipes, which will prevent the proteins from destaining.

STEP	Heated Protocol	Ambient Temperature Protocol
Stain	Rock with 50ml warm (50-65°C) stain 15 min.	Rock with 50ml R250 stain 1 to 2 hours
Destain	Rock with 50ml warm R250 destain with paper towels or lab wipes for 1 to 24 hrs	Rock with 50ml R250 destain 3 to 24 hrs.

3) Staining with “Safe” Stains (Coomassie G250-based stains)

For use with Simply Blue™ and Gel Code™ Blue or Bio-Safe™ stains:

STEP	Heated Protocol	Ambient Temperature Protocol
Fix and Wash	Heat 20% Ethanol at 50°-70°. Rock 5 minutes, 50ml/gel You must do this step.	Rock 30 minutes with 50ml 20% Ethanol per gel
Water Wash	Heat water to circa 80°C. Rock 2 minutes, 100ml/gel. Repeat once. Do not use temperatures higher than 80°C.	Rock 5 minutes with 100ml water per gel. Repeat twice.
Stain	Heat stain to 70°-80°C. Rock 10 minutes, 25ml/gel	Use 25ml stain per gel. Rock 1 hour.
Destain	Heat 10% Ethanol (50ml/gel) at 50°-70°C. Rock 15 minutes - 16 hours with adsorbent paper	Use 100ml 10% ethanol per gel. Rock 2 to 16 hours.
Storage	Add 20ml 20%NaCl in water (w/v) per 100ml of 10% ethanol to preserve results for up to 2 days before drying.	

GENERAL NOTES ON USING COOMASSIE STAINS

One may use either the heated or ambient protocol at any step, but the first step works best when heated as it fixes the proteins better. For the sharpest bands and fastest results, C.B.S. Scientific recommends using the heated protocols. Notes on heating solutions in a microwave:

- 1) Do not heat gel in microwave for temperatures above 70°C, as solution may get too hot and dissolve the agarose.
- 2) For 10% or 20% Ethanol solutions, heat 50ml for 30 seconds on high (600W). The water portion may be heated separately, and alcohol added after heating.
- 3) For water, heat 100ml for 1 minute on high (600W) and then add to the gel.
- 4) Microwave ovens may vary. Check temperature results for your oven and adjust times.

d) Staining with Silver-Stain Solutions


For ammoniacal silver staining protocols (e.g., Pierce Silver SNAP™ or Invitrogen SilverXpress™ Silver Stain kits), fix proteins for 10 minutes with 200ml per gel of 50% methanol/10% acetic acid/20mM sodium bisulfite. For the bisulfite, use Cat. No. LB30001 Antioxidant: 1ml per 200ml fixative. For the SilverXpress kit, use the Tris Glycine protocol using the above fixative, then follow the remaining protocol. For the Silver SNAP kit, use the Tris Glycine protocol using the above fixative, then follow then follow the remaining protocol.

GEL DRYING

Gels may be dried without cracking between cellophane after equilibrating with ClearPAGE™ Drying Solution (HB04510) or with 5% glycerol in water (25ml per gel, rock for 30 minutes). No alcohol is required.

A) Instructions for Gel Drying stained gels using the HB90001 air drying frame:

1. After de-staining, rinse the stained gel once in distilled water (100ml/gel) for 5 minutes.
2. Equilibrate the rinsed gel in ClearPAGE drying solution (HB04510) or water containing 5% glycerol (50ml/gel) for 10 minutes at room temperature.
3. Soak two pieces of cellophane (HB02025) in distilled water.
4. Place the solid frame on 1 or 2 pipette tip boxes.
5. Place one piece of wet cellophane evenly on top, followed by the gel (s), a second piece of wet cellophane, and finally the open frame. (Alternatively, use two open frames for faster gel drying).
6. Carefully smooth away any air bubbles between the cellophane sheets.
7. Clip the frames together, and place on lab bench. Allow to air dry horizontally overnight.
8. When the gel is completely dried, remove the frame and cut away excess cellophane for storage.

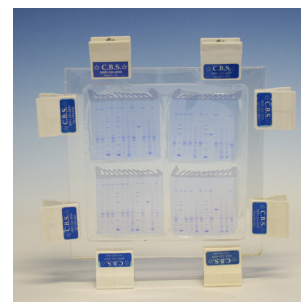

HB90001

B) Instructions for Gel Drying stained gels using the HB90025 air drying frame:

1. After de-staining, rinse the stained gel once in distilled water (100ml/gel) for 5 minutes.
2. Equilibrate the rinsed gel in ClearPAGE drying solution (HB04510) or water containing 5% glycerol (25ml/.gel) for 10 minutes at room temperature, rocking.
3. Soak two pieces of cellophane (HB02520) in distilled water.




4. Place frame #1 onto clamping base, so that frame is flush with clamping platform as shown at right.




5. Place one piece of wet cellophane sheet evenly onto platform.



6. Place up to 4 ClearPAGE™ gels onto cellophane making sure all gels are inside frame perimeter.

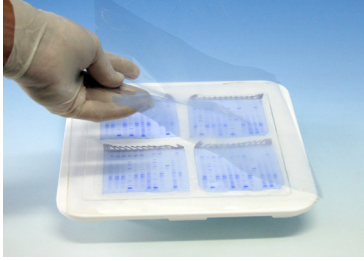


7. Place 2nd damp cellophane sheet on top of gels. Carefully smooth away any air bubbles between sheets.

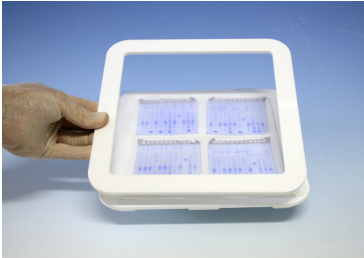


8. Place frame #2 on top of cellophane gel sandwich.



9. Taking advantage of the 8 cut-out areas of the base, clamp frame/ cellophane/gel sandwich ONLY.

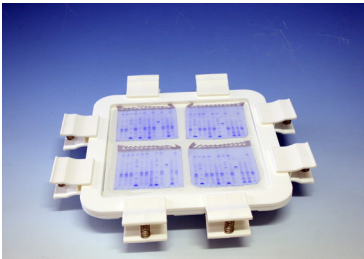

10. Before lifting off clamping base make sure you apply all 8 clamps.

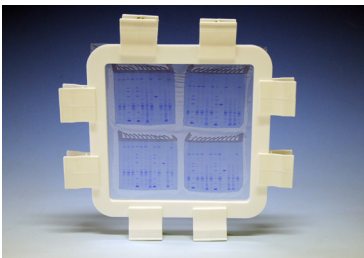

11. Lift off base (as shown) and then place in HORIZONTAL position (not shown) to dry overnight. When gel is completely dried, remove frame and cut away excess cellophane for storage.

GEL BLOTTING


Follow general guidelines as indicated for respective blotting units. As a general recommendation, equilibrate gels (after running) with the diluted transfer buffer for 5 to 10 minutes prior transfer. Clear PAGE™ Transfer Buffer 10x concentrate contains 0.25 M Tris (base), 1.92M Glycine, 1% SDS.

For C.B.S. Scientific's EBU-204, EBX-700, DCX-700, EBU-4000 and EBU-6000 Blotter Buffers:

Component	EBX-700 EBU-204 Tank Blotters	DCX-700 Blotter	EBU-4000 and EBU-6000 Semi-Dry Blotters
ClearPAGE™ Transfer Buffer 10X or 20X dilution (cat. FB82500)	50ml (1:20 dilution)	100ml (1:10 dilution)	10ml (1:10 dilution)
Methanol	200ml	200ml	20ml (for PVDF membrane: reduce 50%)
Ultrapure water	770ml	720ml	72 ml (PVDF: add water for less MeOH)

Typical Blotting conditions	EBX-700, EBU-204 Tank Blotters	DCX-700 or Quadra	EBU-4000 & EBU-6000 Semi-Dry Blotting
Power Supply Setting	50-60V constant	200V constant	25 Volts
Blot time	2-4 hours with stirring. Ex- change cooling blocks after 1.5 hours	1.5 - 2.0 hours with stirring, cooling blocks	30-60 minutes
Expected current	250-300mA	180mA / 1 gel 220mA / 2 gels 440mA / 4 gels	Initial 200-300mA / Final 60- 100mA

Refer to figure below to assemble blotting stack. With cassette wide open assemble components on black side in the following order: buffer saturated sponge pad, gel equilibrated in transfer buffer, buffer saturated transfer membrane, then buffer saturated blotting paper. Smooth with gloved finger or roll with glass rod to be sure no bubbles exist between the gel and the transfer membrane.

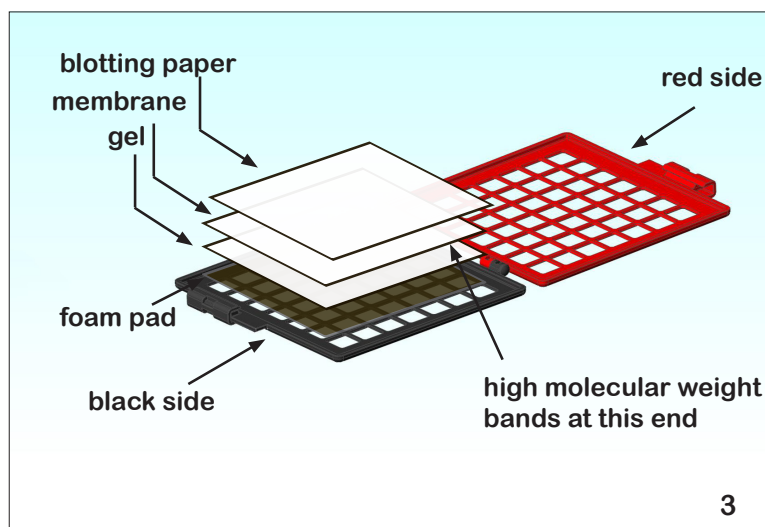

BUFFER FORMULATIONS

As an alternative to the buffers ClearPAGE sells, these formulations may be used to prepare buffers yourself. Use high-quality, low-conductance ingredients. Do NOT use acid or base to adjust the pH!

Standard SDS Running Buffer, 20X for Reduced Samples (FB60500)

Ingredient	MW	Molarity	Qty/Liter
Tricine (free acid)	179.17	0.8M	143.4 g
Tris (free base)	121.14	1.2M	145.2 g
SDS (2%)	288.38	-	20.0 g
Sodium Meta-bisulfite	104.06	50mM	5.0 g
Ultrapure water (fill to)	-	-	1000ml

* pH should be between 8.4 and 8.5 at 25° C.

* For non-reduced samples (especially antibodies), omit the Sodium Meta-bisulfite

Turbo SDS Running Buffer, 20X FB13500

Ingredient	MW	Molarity	Qty/Liter
MPS (free acid)	209.26	0.6M	125.6 g
Tris (free base)	121.14	1.2M	145.2 g
SDS (2%)	288.38	-	20.0 g
Sodium Meta-bisulfite	104.06	50mM	5.0 g
Ultrapure water (fill to)	-	-	1000 ml

* pH should be between 8.3 and 8.4 at 25° C.

* For non-reduced samples (especially antibodies), omit the Sodium Meta-bisulfite

LDS Sample Buffer, 4x -FB31010

Ingredient	MW	Molarity	Qty/Liter
Glycerol (40%)	-	-	400 g
Ficoll-400 (4%)	-	-	40 g
Triethanol amine, pH7.6	149.2	0.8M	120.0
6 N HCL	36.46	-	93.0 g
Lithium Dodecyl Sulfate (4%)	-	-	40 g
EDTA Di-Sodium	372.2	2mM	7.44 g
Brilliant Blue G250 (0.025%)	-	-	0.25g
Phenol Red	-	-	0.25 g
Ultrapure water (fill to)	-	-	1000 ml

* pH should be between 7.7 and 7.8 at 25° C.

**Tris-Glycine-SDS Transfer Buffer, 10X (FB82500) and
ClearPAGE “Classics” Run Buffer, 20X (FB82500)**

Ingredient	MW	Molarity	Qty/Liter
Tris (free base)	121.14	0.25M	30.3 g
SDS (2%)	288.38	-	20.0 g
Glycine	75.07	1.92M	144.1 g
Ultrapure water (fill to)	-	-	1000 ml

* pH should be between 8.4 and 8.6 at 25° C.

Standard DNA/Native Running Buffer, 20X (GB61500)

Ingredient	MW	Molarity	Qty/Liter
Tricine (free acid)	179.17	0.8M	143.4 g
Tris (free base)	121.14	1.2M	145.2 g
Ultrapure water (fill to)	-	-	1000ml

* pH should be between 8.35 and 8.45 at 25° C.

DNA/Native Sample Buffer, 4x (GB33002)

Ingredient	MW	Molarity	Qty/Liter
Glycerol (40%)	-	-	400 g
Ficoll-400 (4%)	-	-	40 g
Triethanol amine, pH7.6	149.2	0.8M	120.0
6 N HCL	36.46	-	93.0 g
EDTA Di-Sodium	372.2	2mM	7.44 g
Brilliant Blue G250 (0.025%)	-	-	0.25g
Phenol Red	-	-	0.25 g
Ultrapure water (fill to)	-	-	1000 ml


* pH should be 7.6 at 25° C.

ClearPAGE TWO-COLOR SDS MARKER

An orange/blue prestained protein marker is also available for SDS PAGE. The recombinant proteins range in size from 7.6 kDa to 195 kDa. The orange markers are 28 kDa and 71 kDa, the other eight bands are blue. ClearPAGE Two-Color SDS Marker is sold ready to use in a 1x LDS sample buffer with 600 μ l per vial (catalog #HM05160). The purity and band sharpness mean the markers can be used for visualizing the separation while running the gel, checking effectiveness on your blot, or as a routine marker in place of other unstained markers.

Recommended volumes to load are:

Application	12-well gels (4.2mm wide)		17-well gels (2.2mm wide) & Marker well on 2D gels	
	load volume	loads/vial	load volume	loads/vial
See during run	10 μ l	60	5 μ l	120
See on blot membrane	4 μ l	150	2 μ l	300
See with stain	4 μ l	150	2 μ l	300

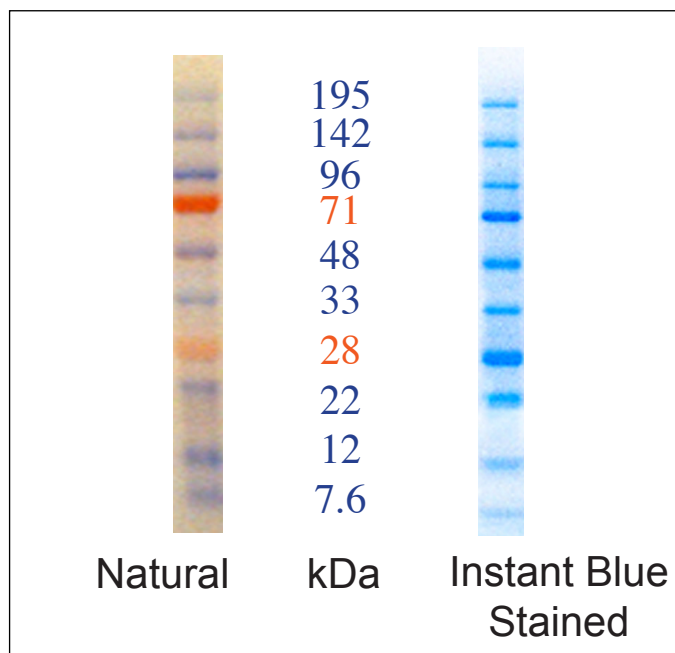

Comparison of Marker on ClearPAGE 12%, 12 well SDS gel with or with Instant Blue stain. Molecular weights are the apparent weights in the ClearPAGE SDS gel system.

Ordering Information for ClearPAGE Precast SDS Gels and Accessories

SPECIAL PROMOTIONS & SAMPLER KITS for SDS PRECAST GELS	
Part number	Description
FM95005	Special offer includes: FREE DCX-700 with purchase of 50 ClearPAGE gels and related buffers
FM91006	Special offer includes: DCX-700 and EPS-300-II Power Supply
FM90002	ClearPAGE Sampler Kit with trial size buffers and 2 SDS gels of customer's choice
FM90005	ClearPAGE SUPER Sampler Kit with trial size buffers and 2 SDS gels of customer's choice, trial size 2-Color SDS marker and trial size Instant Blue Stain
BUFFERS, REDUCING SOLUTIONS, MARKERS & STAINS FOR PROTEIN SEPARATIONS	
Part number	Description
FB31002	ClearPAGE LDS Sample Buffer -2ml
FB31010	ClearPAGE LDS Sample Buffer -10ml
FB32001	ClearPAGE DTT Reducer 10x, 1ml
FB50053	ClearPAGE SDS-Non-Reducing Running Buffer - 50ml, 20x
FB50500	ClearPAGE SDS-Non-Reducing Running Buffer - 500ml, 20x
FB60053	ClearPAGE SDS-Reducing Running Buffer -50ml, 20x
FB60500	ClearPAGE SDS-Reducing Running Buffer - 500ml, 20x
FB14053	ClearPAGE TURBO Reducing Running Buffer- 50ml, 20x
FB14500	ClearPAGE TURBO Reducing Running Buffer- 500ml, 20x
FB82500	ClearPAGE TGS Transfer Buffer-500ml
HG73010	ClearPAGE Instant Blue Stain, 1 liter, enough for 40-50 mini-gels
HG73050	ClearPAGE Instant Blue Stain, 50ml, trial size
HM05160	2-Color SDS Marker, 0.6ml, recombinant protein range: 7.6kDa - 195kDa
PRECAST GELS-SDS- 10cm x 10cm for use with TEO Buffer (tris-trycine with tri-ethanol amine)	
Part number	Description
FK00812-10	ClearPAGE SDS Gel 8%, 12 well, package of 10
FK00827-10	ClearPAGE SDS Gel 8%, 17 well, package of 10
FK01001-10	ClearPAGE SDS Gel 10%, 1 well, package of 10
FK01012-10	ClearPAGE SDS Gel 10%, 12 well, package of 10
FK01027-10	ClearPAGE SDS Gel 10%, 17 well, package of 10
FK01201-10	ClearPAGE SDS Gel 12%, 1 well, package of 10
FK01212-10	ClearPAGE SDS Gel 12%, 12 well, package of 10
FK01227-10	ClearPAGE SDS Gel 12%, 17 well, package of 10
FK01612-10	ClearPAGE SDS Gel 16%, 12 well, package of 10
FK01627-10	ClearPAGE SDS Gel 16%, 17 well, package of 10
FK40812-10	ClearPAGE SDS Gel 4-8%, 12 well, package of 10
FK40827-10	ClearPAGE SDS Gel 4-8%, 17 well, package of 10
FK41212-10	ClearPAGE SDS Gel 4-12%, 12 well, package of 10
FK41227-10	ClearPAGE SDS Gel 4-12%, 17 well, package of 10
FK42001-10	ClearPAGE SDS Gel 4-20%, 1 well, package of 10
FK42002-10	ClearPAGE SDS Gel 4-20%, 2D, package of 10
FK42012-10	ClearPAGE SDS Gel 4-20%, 12 well, package of 10
FK42027-10	ClearPAGE SDS Gel 4-20%, 17 well, package of 10
FK82001-10	ClearPAGE SDS Gel 10-20%, 1 well, package of 10
FK82012-10	ClearPAGE SDS Gel 10-20%, 12 well, package of 10
FK82017-10	ClearPAGE SDS Gel 10-20%, 17 well, package of 10
PRECAST GELS- SDS- 10cm x 8cm- for use with TEO Buffer (tris-trycine with tri-ethanol amine)	
Part number	Description
BRGK-100	Replacement Gasket adapts Bio Rad's Mini-Protean II for use with 8cm x10cm ClearPAGE gels
BK01212-10	ClearPAGE SDS Gel 12%, 12 well, package of 10
BK01227-10	ClearPAGE SDS Gel 12%, 17 well, package of 10
BK40812-10	ClearPAGE SDS Gel 4-8%, 12 well, package of 10
BK40827-10	ClearPAGE SDS Gel 4-8%, 17 well, package of 10
BK42012-10	ClearPAGE SDS Gel 4-20%, 12 well, package of 10
BK42027-10	ClearPAGE SDS Gel 4-20%, 17 well, package of 10
BLOTTING & DRYING ACCESSORIES for SDS GELS	
Part number	Description
FB82500	ClearPAGE TGS Transfer Buffer - 500ml
HB90025	4-place Drying Frame with "cut-out", 25cm x 25cm with 8 clamps, 20 pieces of cellophane, and 125mls of drying solution
HB90013	4-place Drying Frame, with 8 clamps and 20 pieces of cellophane
HB01320	Precut cellophane - 13 x 13 cm, package of 20 sheets
HP19001	ClearPAGE Blot Sandwich NC 90x85mm 2/pk. 0.2µ Nitrocellulose/1mm blotting paper
HP19020	ClearPAGE Blot Sandwich NC 90x85mm 20/pk. 0.2µ Nitrocellulose/1mm blotting paper
HP29301	ClearPAGE Blot Sandwich PVDF 90x85mm 2/pk. 0.2µ PVDF/1mm blotting paper
HP29320	ClearPAGE Blot Sandwich PVDF 90x85mm 20/pk. 0.2µ PVDF/1mm blotting paper
HB04510	ClearPAGE Gel Drying Solution - 1 liter
HB02520	Precut Cellophane - package of 20 sheets
HB90001	Drying frame with 8 clamps and 20 sheets of cellophane

Ordering Information for ClearPAGE DNA/Native Precast Gels & Accessories

SPECIAL PROMOTIONS & SAMPLER KITS for DNA/NATIVE PRECAST GELS	
Part number	Description
FM90003	ClearPAGE Sampler Kit with trial size buffers and 2 DNA/Native gels of customer's choice
BUFFER SOLUTIONS FOR DNA/NATIVE PROTEIN SEPARATIONS	
Part number	Description
FB82500	ClearPAGE TGS Transfer & Running Buffer - 500ml
GB33010	ClearPAGE DNA/Native Buffer -10ml
GB61053	ClearPAGE DNA/Native Running Buffer-50ml
GB61500	ClearPAGE DNA/Native Running Buffer- 0.5 liters
PRECAST GELS-DNA/NATIVE	
Part number	Description
GN01012-10	ClearPAGE DNA/Native Gel 10%, 12 well, package of 10
GN01027-10	ClearPAGE DNA/Native Gel 10%, 17 well, package of 10
GN20812-10	ClearPAGE DNA/Native Gel 2-8%, 12 well, package of 10
GN20827-10	ClearPAGE DNA/Native Gel 2-8%, 17 well, package of 10
GN32012-10	ClearPAGE DNA/Native Gel 3-20%, 12 well, package of 10
GN32027-10	ClearPAGE DNA/Native Gel 3-20%, 17 well, package of 10
DRYING ACCESSORIES for DNA/NATIVE GELS	
Part number	Description
HB04510	ClearPAGE Gel Drying Solution - 1 liter
HB02520	Precut Cellophane - package of 20 sheets
HB90001	Drying frame with 8 clamps
HB90025	4-place Drying Frame with "cut-out", 25cm x 25cm with 8 clamps, 20 pieces of cellophane, and 125mls of drying solution

Ordering Information for ClearPAGE "Classics" TGS Precast Gels & Accessories

SPECIAL PROMOTIONS & SAMPLER KITS for ClearPAGE "Classics" TGS PRECAST GELS	
Part number	Description
FM90003	ClearPAGE Sampler Kit with trial size buffers and 2 DNA/Native gels of customer's choice
BUFFER SOLUTIONS for TGS gels	
Part number	Description
FB82500	ClearPAGE TGS Transfer and TGS ClearPAGE Classics Running Buffer - 500ml
PRECAST GELS-DNA/NATIVE	
Part number	Description
TG01212-10	ClearPAGE "Classics" TGS Gel 12%, 12 well, package of 10
TG01217-10	ClearPAGE "Classics" TGS Gel 12%, 17 well, package of 10
TG42012-10	ClearPAGE "Classics" TGS Gel 4-20%, 12 well, package of 10
TG42017-10	ClearPAGE "Classics" TGS Gel 4-20%, 17 well, package of 10
TG81612-10	ClearPAGE "Classics" TGS Gel 8-16%, 12 well, package of 10
TG81617-10	ClearPAGE "Classics" TGS Gel 8-16%, 17 well, package of 10
DRYING ACCESSORIES for ClearPAGE Classics	
Part number	Description
HB04510	ClearPAGE Gel Drying Solution - 1 liter
HB02520	Precut Cellophane - package of 20 sheets
HB90001	Drying frame with 8 clamps
HB90025	4-place Drying Frame with "cut-out", 25cm x 25cm with 8 clamps, 20 pieces of cellophane, and 125mls of drying solution

Note: ClearPAGE "Classics" TGS Precast gels have a separate set of instructions

Ordering Information for Related Products

DUAL COOL SYSTEM	
Part number	Description
DCX-700	Dual Cool System, CE. 2 place Combo System. Fits precast gels or glass plate dimensions of 10x10cm or 10x8cm(h). Kit includes: lower reservoir, safety cover with attached leads, core, 2 blotting cassettes with sponge pads, 2 freezer blocks, single gel adapter plate and instruction manual
QNX-700	Quadra System, CE, 4 place Combo System without cooling. Two cores accomodate 4 gels or 4 blotting cassettes. Fits 4 precast gels or 4 sets of glass plates with dimensions of 10x10cm or 10x8cm(h). Kit includes: lower reservoir, safety cover with attached leads, 2 cores, 4 blotting cassettes with sponge pads, and instruction manual
QNC-700	Quadra System, CE, 4 place Combo System with double-sided cooling labyrinth. Two cores accomodate 4 gels or 4 blotting cassettes. Fits 4 precast gels or 4 sets of glass plates with dimensions of 10x10cm or 10x8cm(h). Kit includes: lower reservoir, safety cover with attached leads, 2 cores, 4 blotting cassettes with sponge pads, and instruction manual
4-PLACE BLOTTING SYSTEM	
Part number	Description
EBX-700	4-Place Blotter, CE. Kit includes: lower reservoir, safety cover with attached leads, 4 blotting cassettes, 4 sponge pads, 4 freezer blocks, anode, cathode, and instruction manual
POWER SUPPLIES	
Part number	Description
EPS-300-II	Mini Power Supply with 24hr timer, CV or CC, 10-300V, 110V/60Hz, current range: 4-500mA, 90 watts
EPS-300-IV	Mini Power Supply with 24hr timer, CV or CC, 10-300V, 220V/50Hz, current range: 4-500mA, 90 watts
EPS-200-II	High Current Electroblotting Power Supply with 24hr timer, CV or CC, 5-200V, 110V/50-60Hz, current range: 4-2000mA, 200 watts
EPS-200-IV	High Current Electroblotting Power Supply with 24hr timer, CV or CC, 5-200V, 220V/50-60Hz, current range: 4-2000mA, 200 watts

ClearPAGE™ are trademarks of C.B.S. Scientific Company
Pierce™, Silver Snap™, and Gel Code™ are trademarks of Pierce Biotechnology, Inc.
Invitrogen™, Simply Blue™, SilverXpress™, XCell II™, and XCellSureLock™ are trademarks of Invitrogen Corp.
Coomassie™ is a trademark of Imperial Chemical Industries PLC.

www.cbsscientific.com
phone: 858-755-4959
tollfree: 800-243-4959
Fax: 858-755-0733

C.B.S. SCIENTIFIC
COMPANY, INC

www.cbsscientific.com

phone: 858-755-4959

tollfree: 800-243-4959

Fax: 858-755-0733

C.B.S.★SCIENTIFIC
C O M P A N Y , I N C