
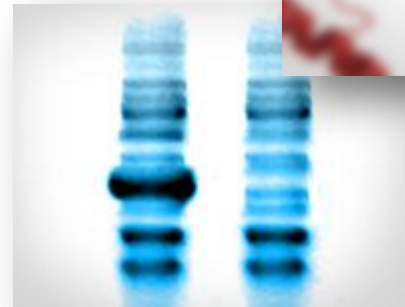
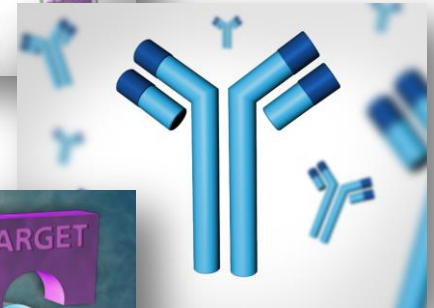


POROS[®] CEX Slides: HS and XS

April 2014

Overview

- CEX introduction: POROS® HS and XS
 - **There are lots of slides that have duplicate information. You will need to chose the ones you would like to present.**
- Applications and Characterization Data
 - Applications
 - Viral clearance
 - Pressure vs flow curves
 - Particle size analysis
- How to optimize
- Case Studies:
 - POROS® HS process redesign
 - POROS® HS in HTP screening and scale up study
 - Rethinking polish applications: CEX in flow through/overload mode
- Summary

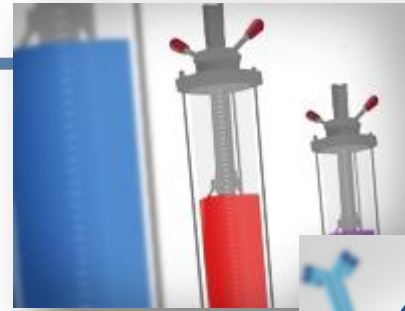

POROS® Process Chromatography Products



- **Cation Exchange**

- HS 50 Strong CEX
- XS Strong CEX

- **Anion Exchange**

- HQ 50 Strong AEX
- PI 50 Weak AEX
- D 50 Weak AEX


- **Affinity**

- MabCapture A rProtein A
- A 50 rProtein A (legacy)
- Heparin Heparin
- OH50 Activated Affinity



POROS[®] HS delivers optimal dynamic binding capacity with high resolution




High dynamic binding capacity over a wide range of conditions

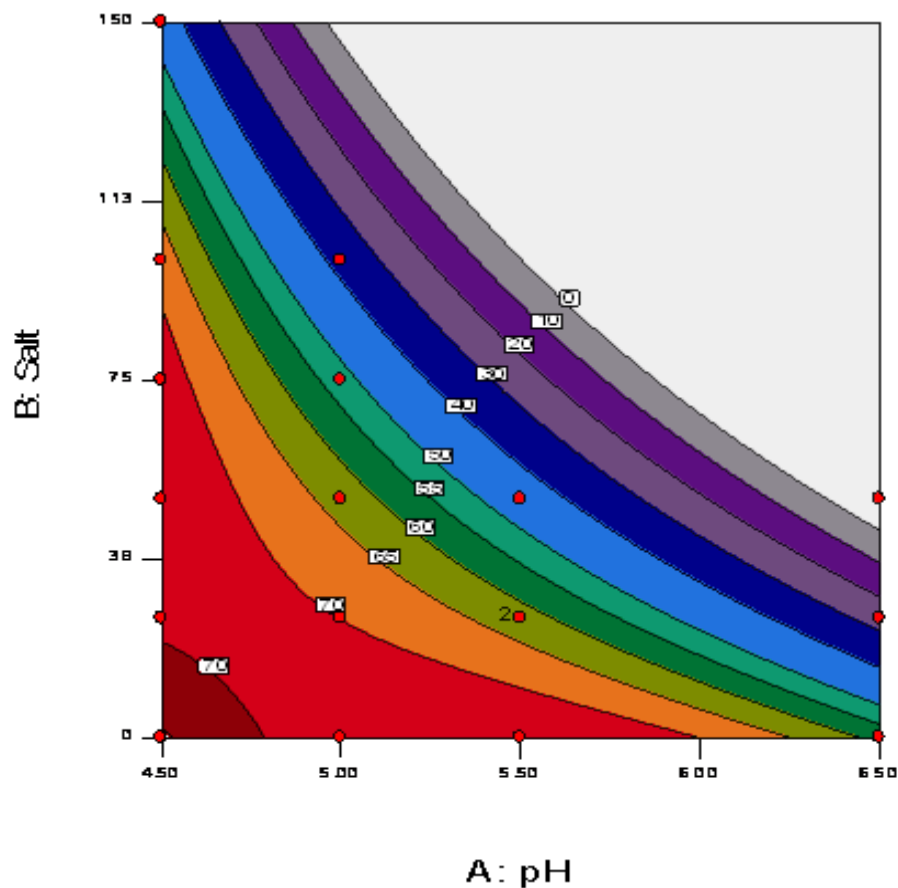

High resolution maintained as linear flow rate increases

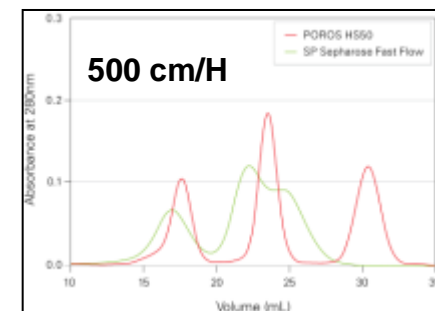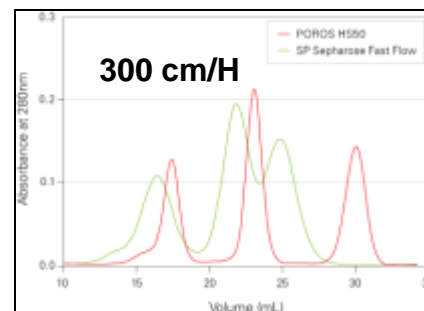

IgG Capacity at 5% Breakthrough vs Linear Flowrate

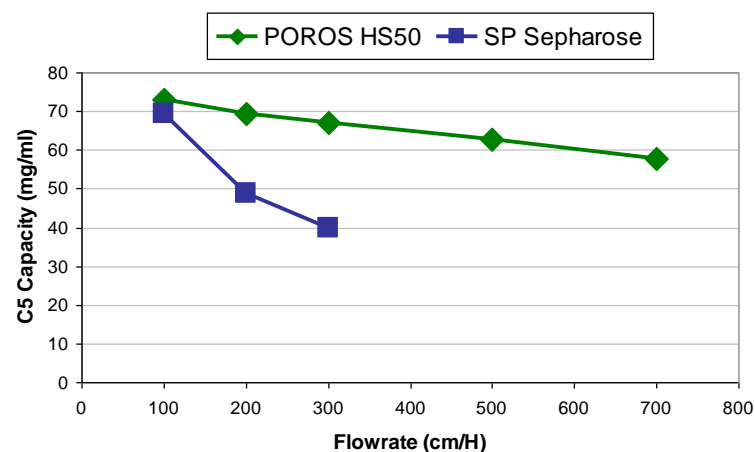


POROS[®] XS delivers optimal dynamic binding capacity with high resolution





High capacity over a broad range of process conditions

IgG Binding Capacity vs. Salt Concentration and pH
Capacity at 5% Breakthrough

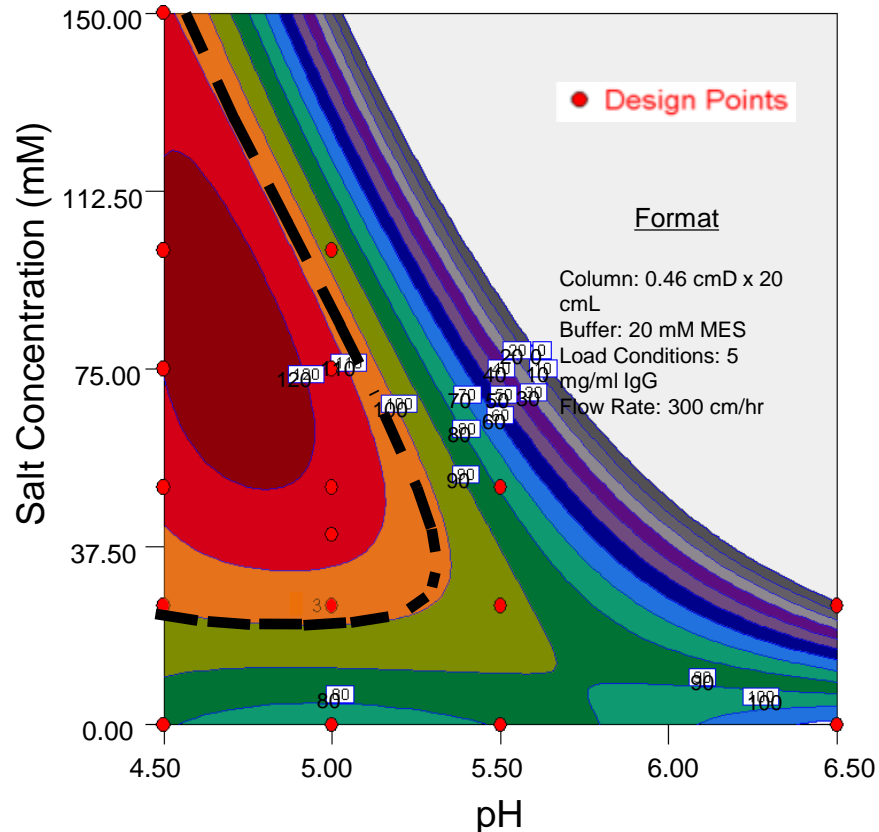

High resolution maintained as linear flow rate increases

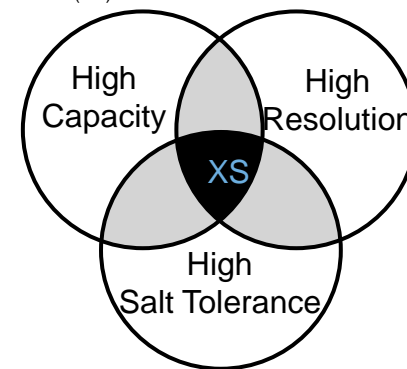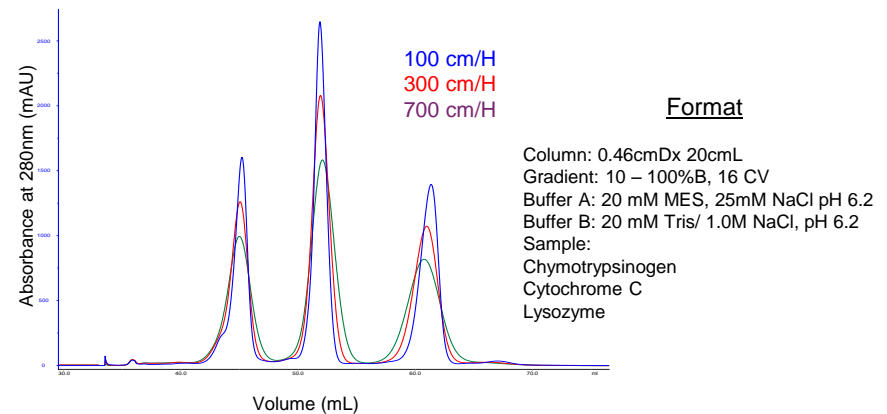


POROS[®] XS delivers optimal dynamic binding capacity with high resolution





High capacity over a broad range of process conditions

IgG Binding Capacity vs. Salt Concentration and pH
Capacity at 5% Breakthrough

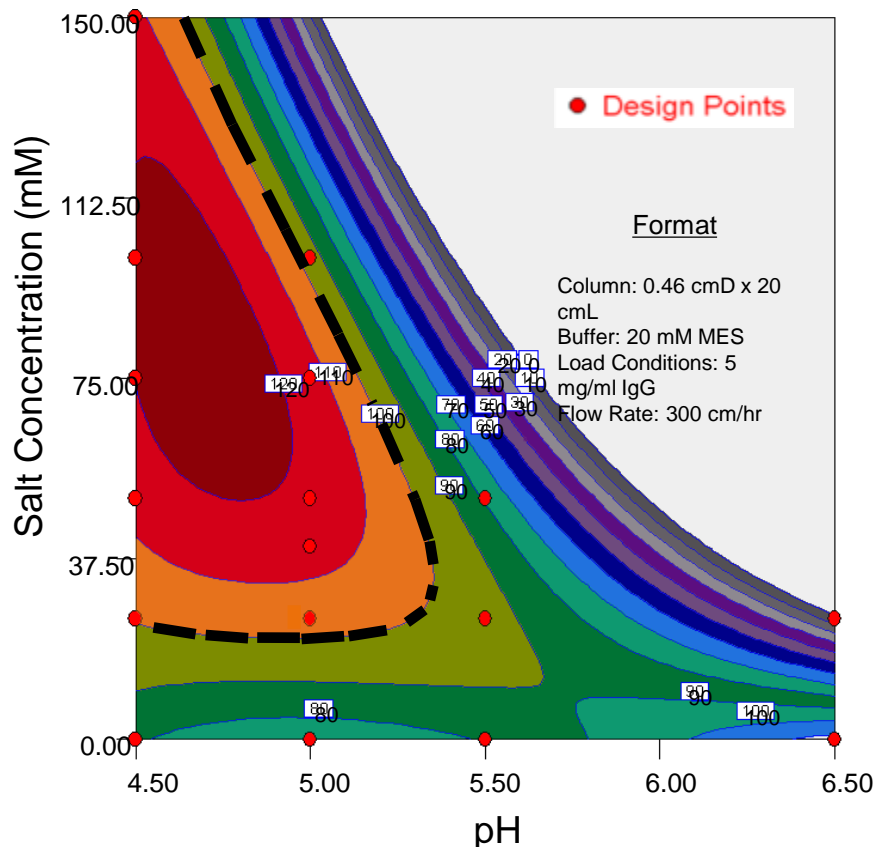

High resolution maintained as linear flow rate increases


Capacity of Cation Exchange Resins:
100 mM NaCl, pH 4.5, 300 cm/H

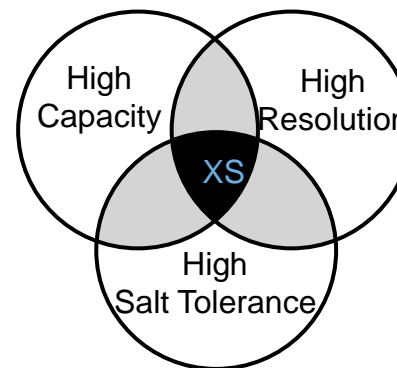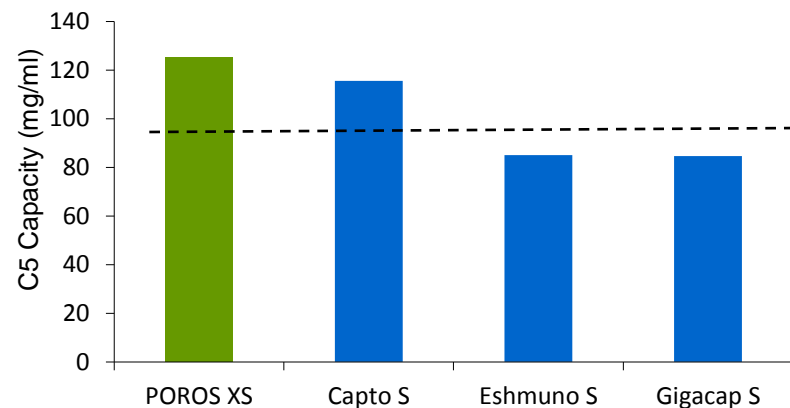



Superior Resolution Independent of Flow Rate

POROS[®] resins maintain resolution as flow rate increases

Separation of Lysozyme, Chymotrypsinogen and Cytochrome C
Column Format: 4.6mmD x 100mmL
Elution Gradient: 15 Column Volumes 0 to 0.5M Sodium Chloride, 20mM MES, pH 6.2

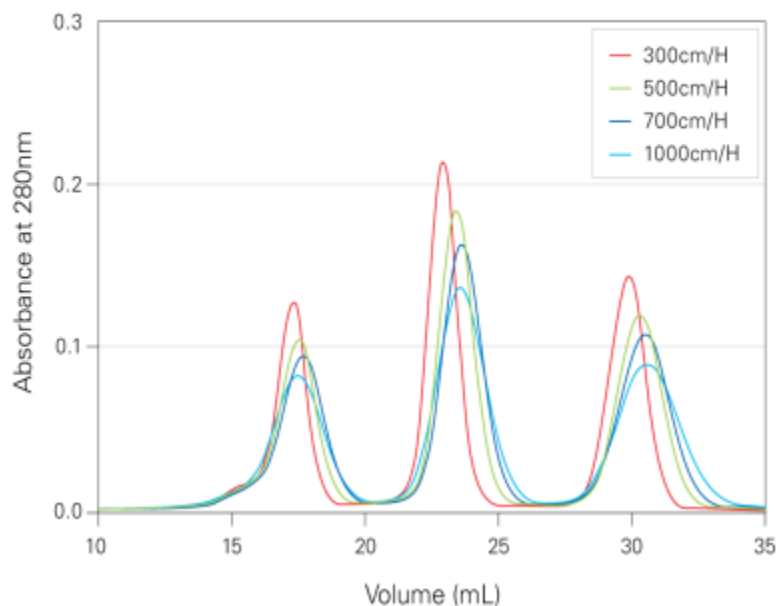

Conventional resins lose resolution as flow rate increases

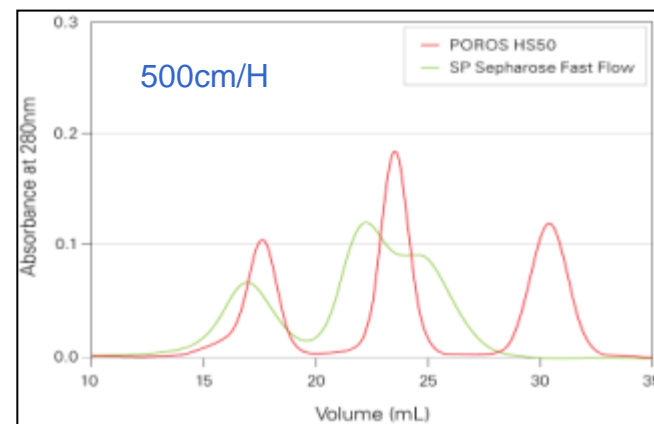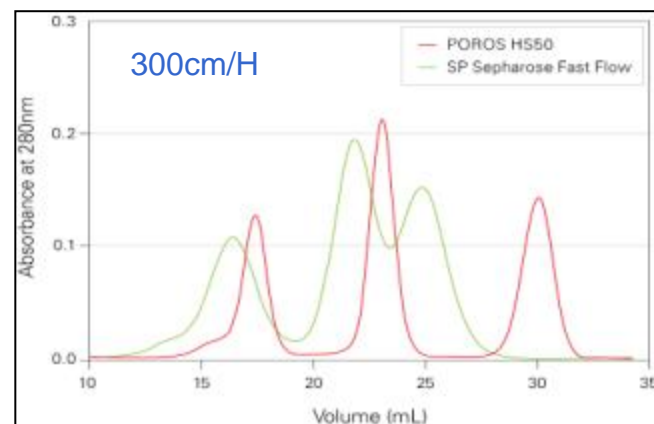

Ion Exchange Chromatography

The workhorse of the purification process

IEX Resin Type	Mode of Operation	Key Applications
Cation Exchange	Bind/elute	Polish of many biomolecules (Mabs, VLP/viruses, fusion proteins, high pI rProteins)
	Overload/ Flow through	Polish for Mabs by binding impurities under normal B/E conditions: impurity removal (aggregates, HCP, DNA, viruses)
Anion Exchange	Bind/elute	Polish of many biomolecules (low pI rProteins, DNA, viruses, plasmids)
	Flow through	Polish for Mabs: binds impurities (DNA, viruses, HCP, aggregates, endotoxin)

POROS® XS

A High Performance Strong Cation Exchanger

- Design Goals:
 - Increase dynamic binding capacity of POROS HS
 - Retain all the beneficial performance attributes of POROS HS
 - resolution, pressure v flow, chemical stability and physical stability
 - Co-developed with industry based biopharma
 - Target flow rate 100-400 cm/hr
 - Increased salt tolerance
- Optimized base bead characteristics
 - Increased total pore volume
 - Decrease average pore size
 - Increased ligand density



POROS[®] XS vs HS

- The material of construction for the POROS XS and HS base beads is polystyrene-divinylbenzene. The sulphopropyl surface charge for both products is derived from the same negatively charged polymer. In addition, POROS XS and HS both have a 50um average particle size and very similar particle size distributions.
- In the development of POROS XS the based bead morphology and surface functionalization was optimized to drive higher IgG dynamic binding capacity under higher conductivity conditions (salt tolerance). In relation to HS, the average pore size for XS is lower (HS = 1600A and XS = 1100A) and the pore volume is higher (HS = x and XS = x + 18%). Using DOE, the sulphopropyl functionalization process was optimized to drive higher ligand density onto the bead surface and the pore morphology was optimized to yield the highest accessibility to the surface charge. Ligand density is assessed with an ionic capacity measurement. Having high accessible ionic capacity leads to high dynamic binding capacity and improved salt tolerance. The ionic capacity for POROS XS is 40-50% higher than the ionic capacity for POROS HS.



POROS® XS Chromatography Resin Overview




- Increased Process Flexibility and Throughput
- Improved Yield
- Smaller Column Sizes
- Linear and Predictable Scalability
- Easier Handling and Packing

Customer Requirements	XS Delivers
High capacity	>100 mg/mL 5% breakthrough capacity for monoclonal antibodies and recombinant proteins under a wide range of process conditions
High resolution	Optimized 50 µm particle size for improved impurity clearance and better yield
High salt-tolerant protein capacity	>100 mg/mL dynamic binding capacity at 15 mS/cm NaCl concentration
Low back pressure	2.0 bar at 800 cm/H in 20 cm length column
Rigid polymeric bead	Incompressible bead with robust physical stability
A high performance CEX resin that combines all of the above	New industry standard – major improvement on industry proven chromatography platform that provides:

POROS[®] XS Chromatography Resin

Design space for optimal dynamic binding capacity

IgG Binding Capacity vs. Salt Concentration and pH
Capacity at 5% Breakthrough

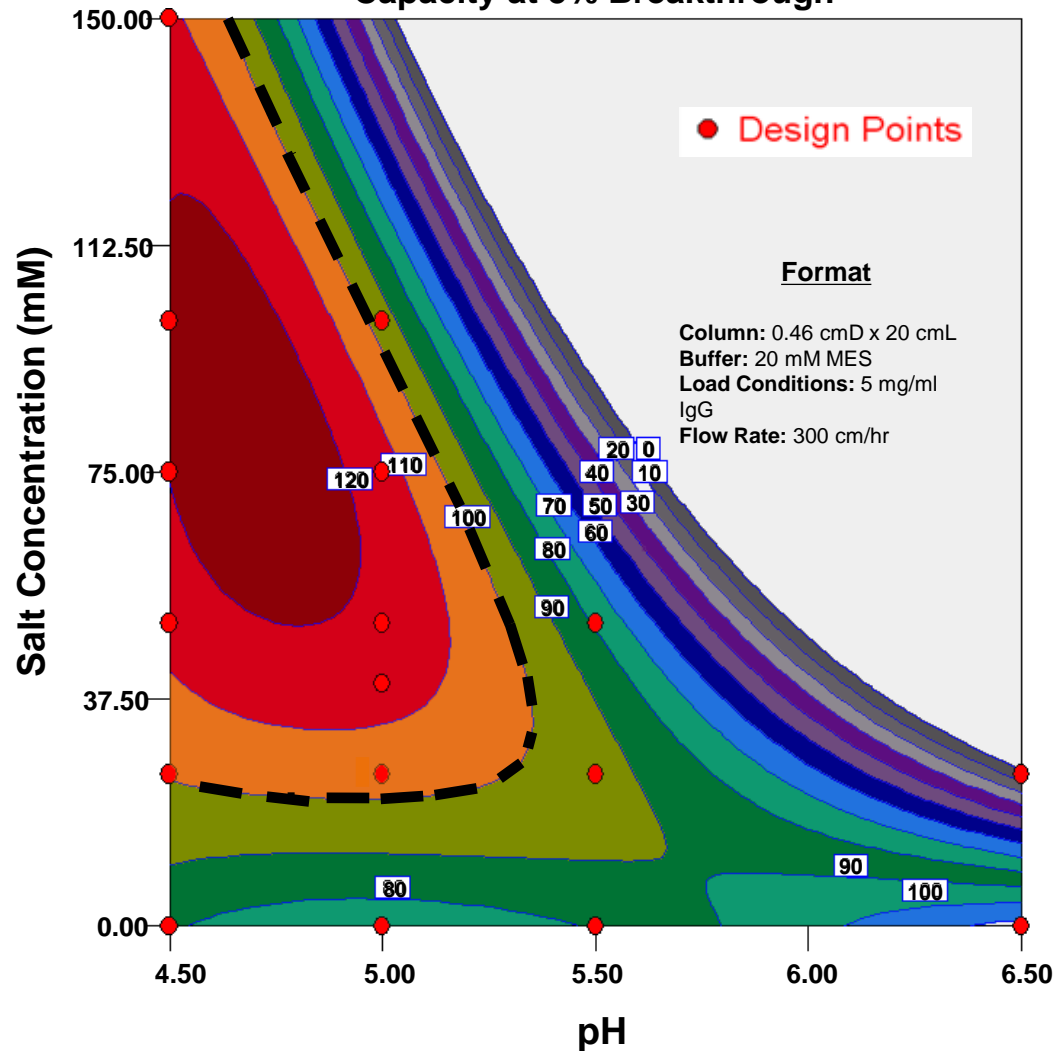


POROS[®] XS
delivers high
capacity over a
broad range of
process
conditions

POROS[®] XS Chromatography Resin

Protein Capacity


IgG Binding Capacity vs. Salt Concentration

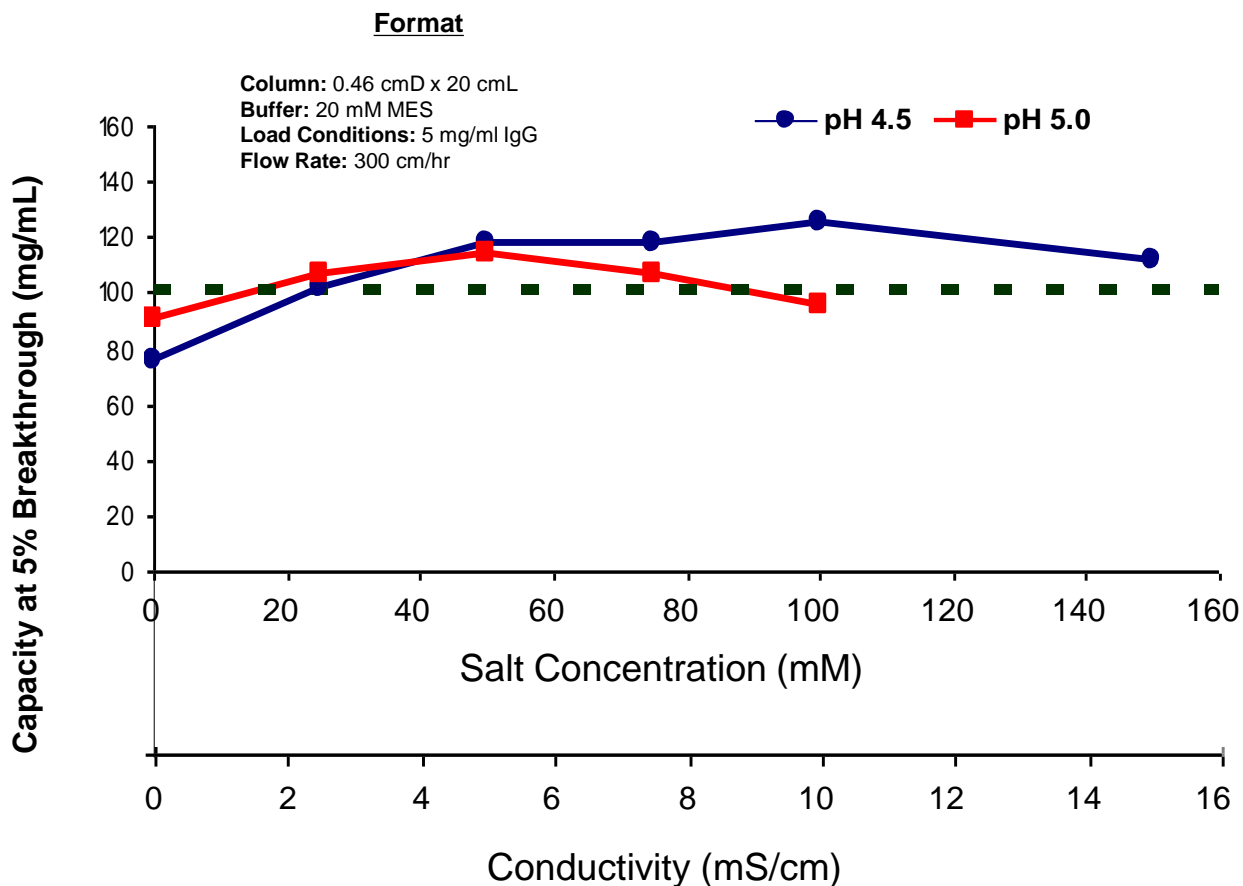

- Achieved greater than 100 mg/ml capacity over a broad range of salt concentrations
- Increased salt tolerance for improved process flexibility

POROS[®] XS IgG Binding Capacity with Increased Flow Rate

POROS XS: IgG Capacity vs Flowrate

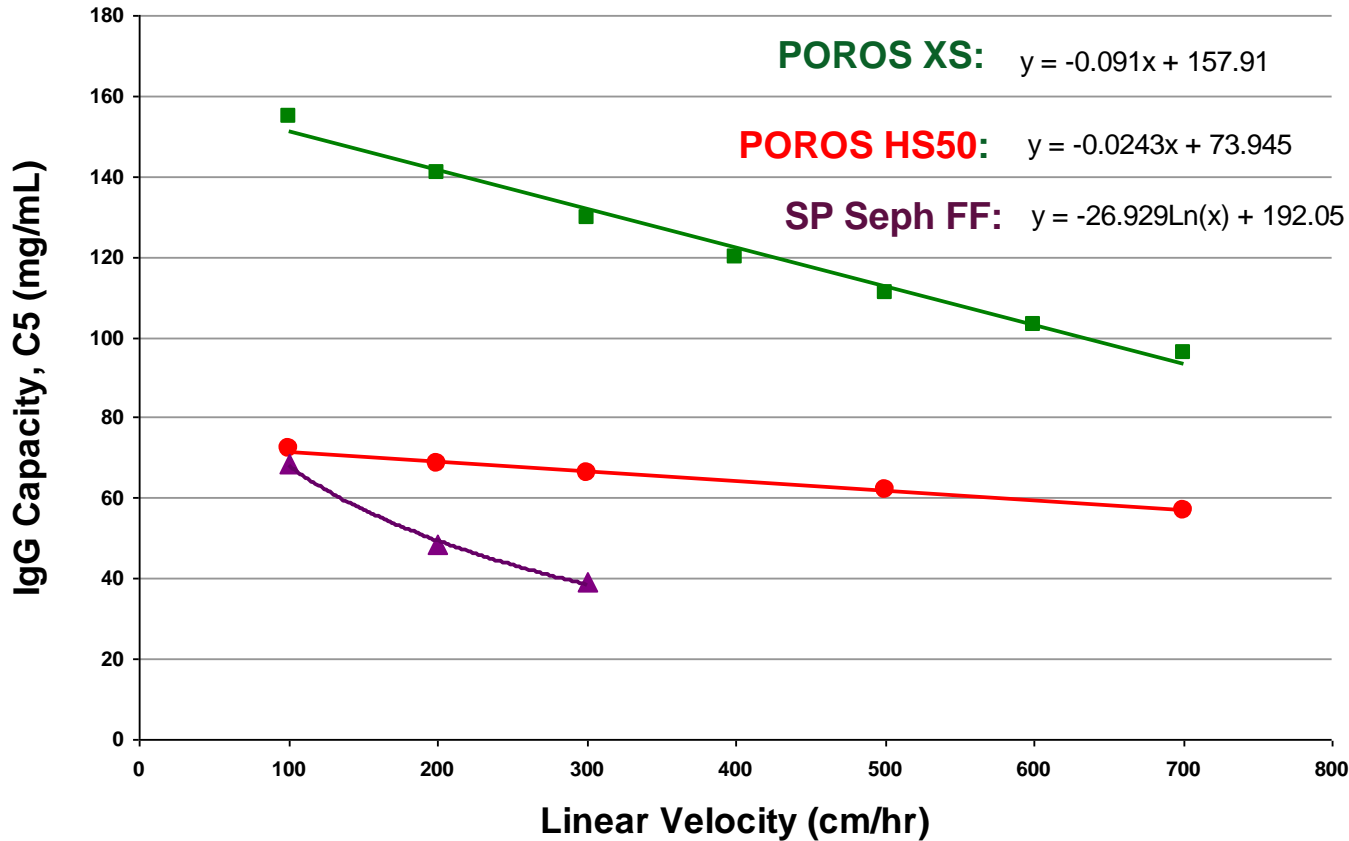

Format

Column: 0.46cmDx 20cmL

Buffer: 20mM MES, 40mM NaCl
pH 5.0

Load Condition: 5 mg/ml
Polyclonal IgG in Equil

Flow Rate: 300 cm/hr

Obtain high
linear binding
capacity over a
large flow rate
range

life
technologies

POROS[®] XS Lysozyme Binding Capacity

POROS[®] XS has >100 mg/ml binding with different proteins with increased salt

Format


Column: 0.46cmDx 20cmL

Buffer: 20mM MES pH 6.2

Load Condition: 7-8 mg/ml
Lysozyme in Equil


Flow Rate: 300 cm/hr

**Lysozyme Dynamic Binding Capacity with
No Salt and 40mM NaCl at 5 and 50% Breakthrough**

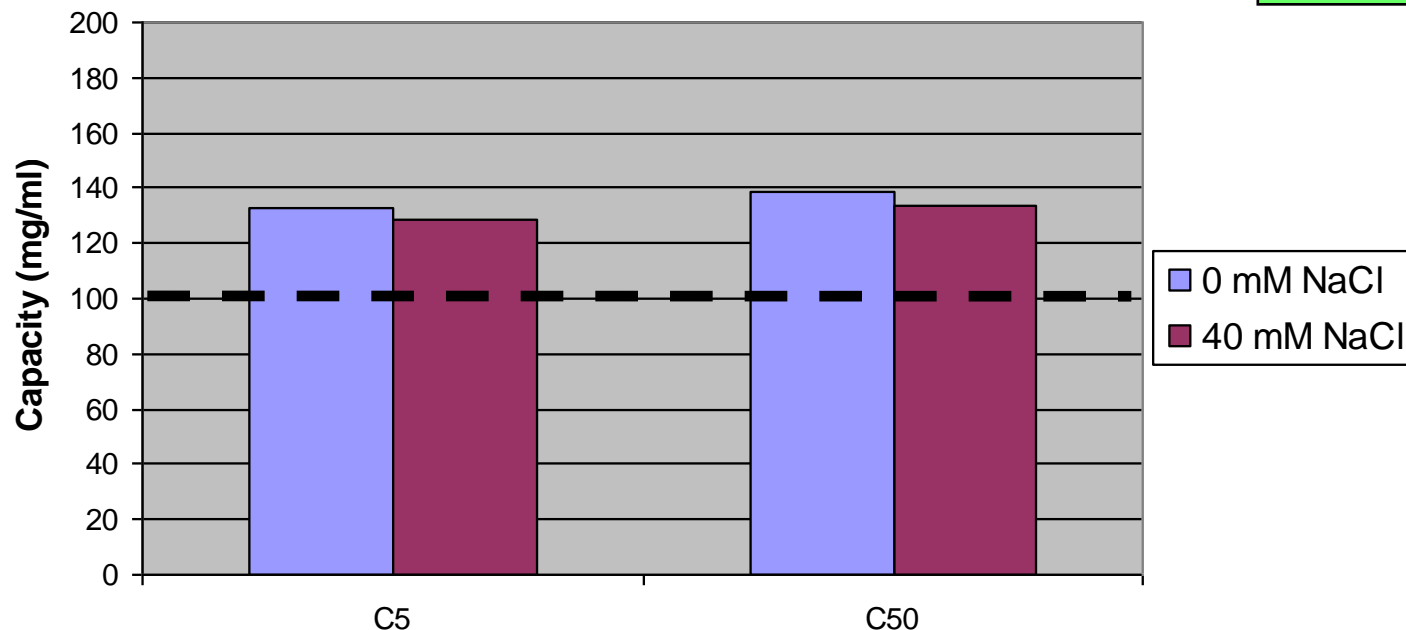

life
technologies

POROS[®] XS Lysozyme Binding Capacity with Increased Flow Rate

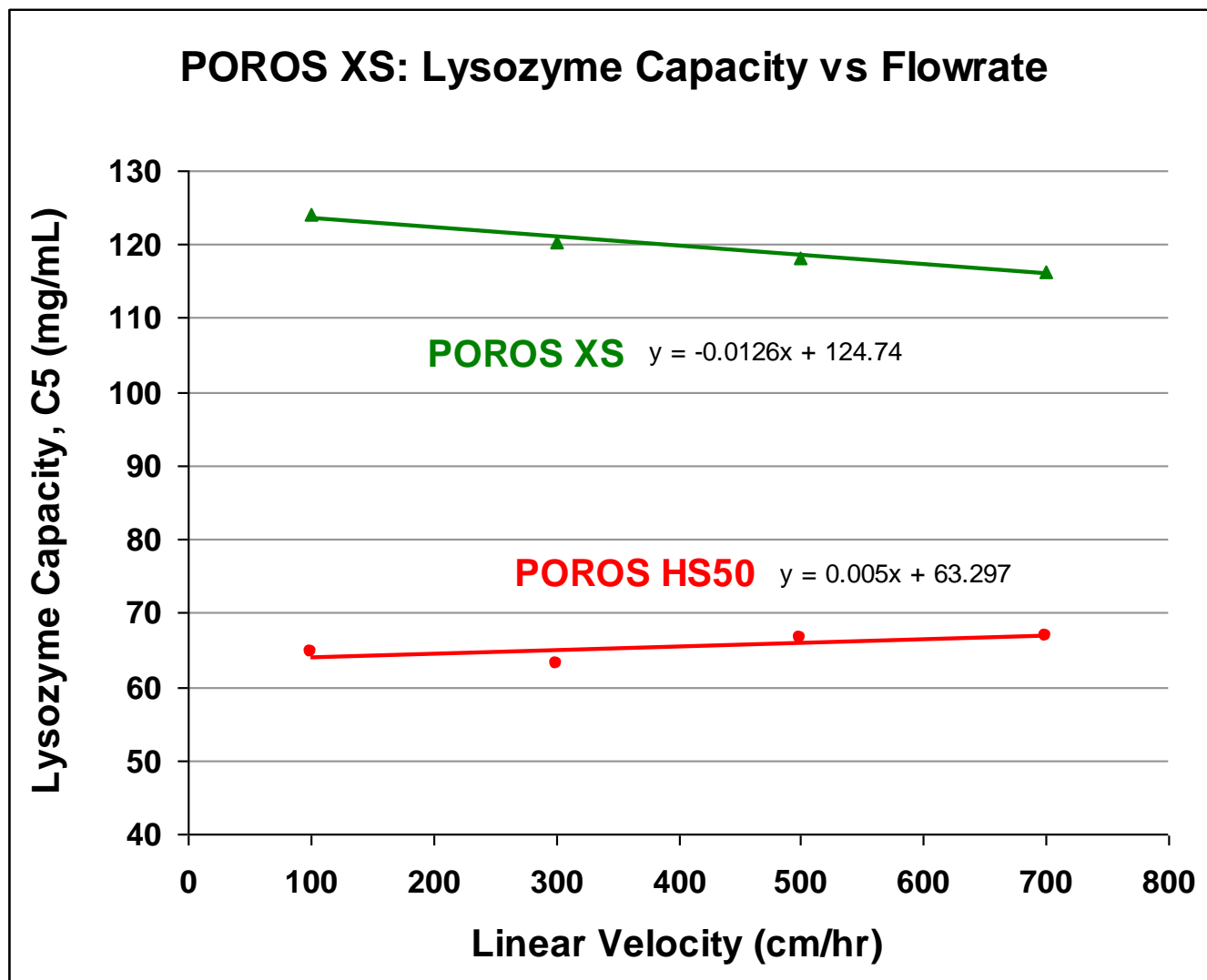

Format

Column: 0.46cmDx 20cmL

Buffer: 20mM MES pH 6.2

Conductivity: 1.1 ms/cm

Load Condition: 3 mg/ml
Lysozyme in Equil


Flow Rate: variable

Obtain high
binding
capacity over
a large flow
rate range

life
technologies

POROS[®] Chromatography Resin

Product Attributes

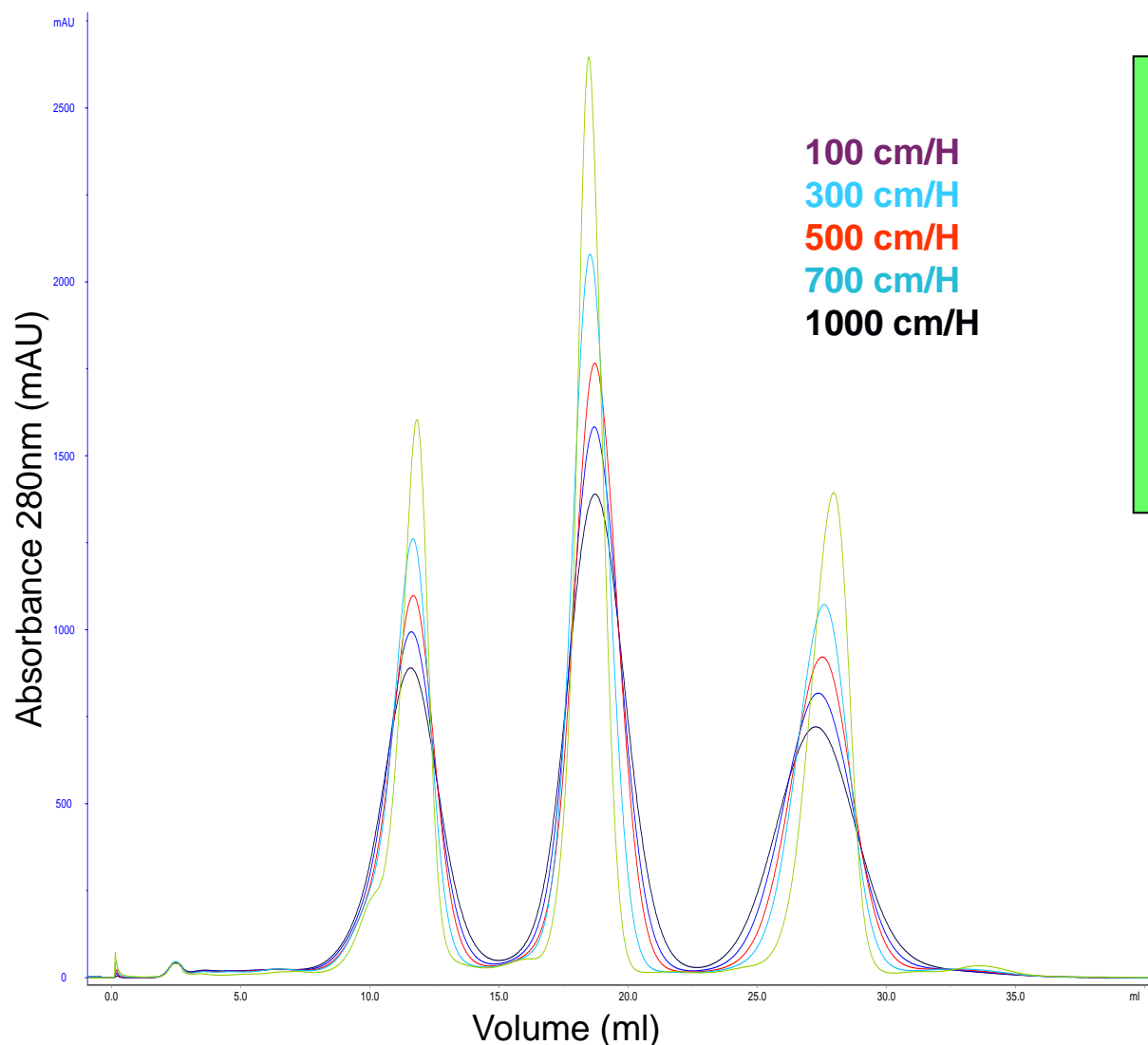

Format

Column: 0.46cmDx 20cmL, 3.32 ml

Gradient: 10 – 100%B, 16 CV

Buffer A: 20 mM MES, 25mM NaCl pH 6.2

Buffer B: 20 mM MES/ 1.0M NaCl, pH 6.2

Sample:

Chymotrypsinogen A 5mg/mL, pI 9.4

Cytochrome C 8mg/mL, pI 10.4

Lysozyme 4mg/ml, pI 11.3

Total Protein Loaded: 5.175mg, 1.56mg/ml

Resolution is
maintained on
POROS[®] resin as
linear flow rate
increases

life
technologies

POROS® XS Characterization

Protein separation compared to POROS® HS

Separation Conditions

Column = 1cmD x 20cmL

Buffer A = 20mM MES, 25mM NaCl pH 6.2


Buffer B = 20mM MES, 1M NaCl pH 6.2

Elution = 10% B – 50% B 7.5CV

Flow Rate = 300cm/H

Protein Mix =

Chymotrypsinogen A
Cytochrome C
Lysozyme

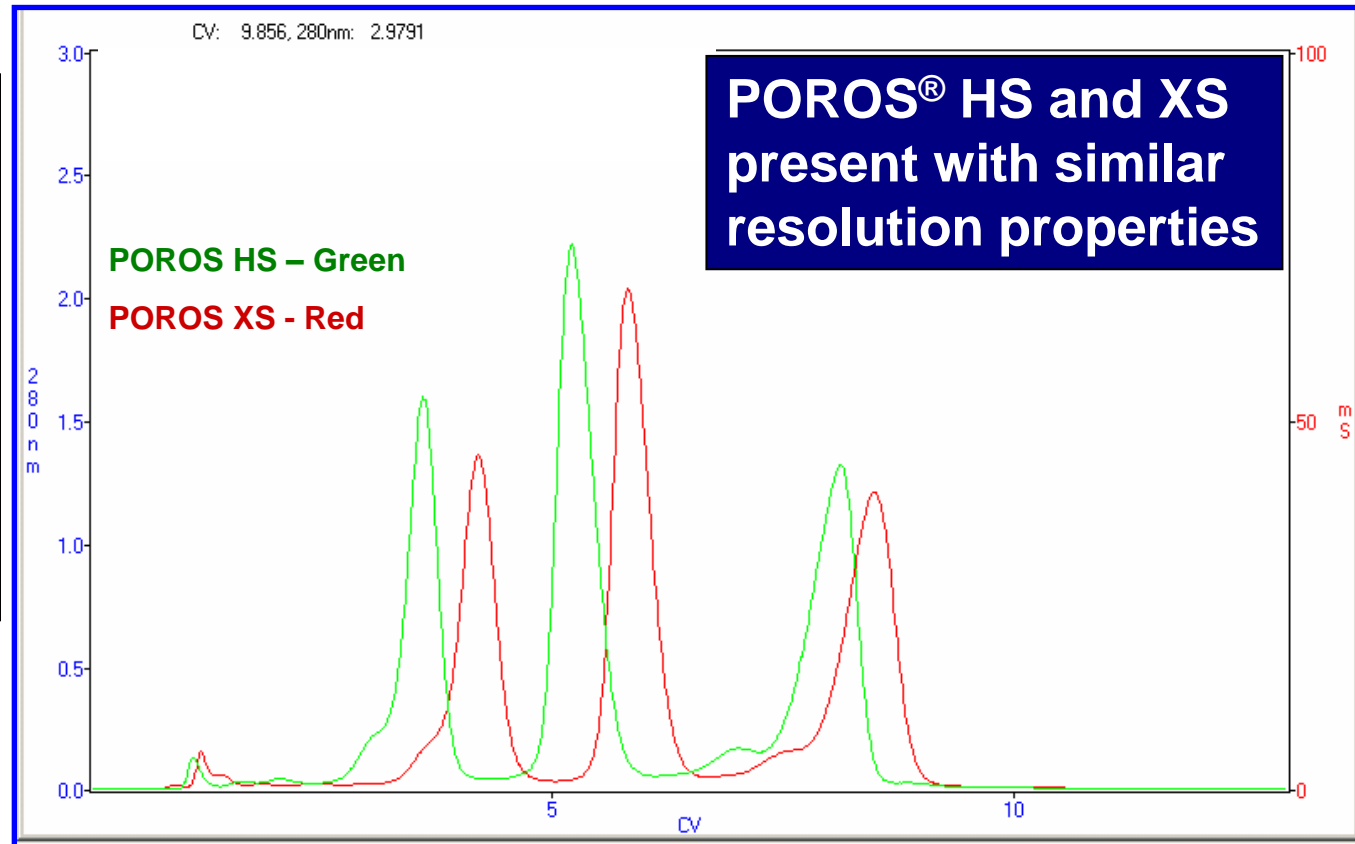

POROS[®] XS Characterization

Protein separation across the Ionic Capacity range

Separation Conditions

Column = 0.46cmD x 20cmL

Buffer A = 20mM MES, 25mM NaCl pH 6.2

Buffer B = 20mM MES, 1M NaCl pH 6.2

Elution = 10% B – 50% B 7.5CV

Flow Rate = 300cm/H

Protein Mix =

Chymotrypsinogen A 5mg/mL, pI 9.4

Cytochrome C 8mg/mL, pI 10.4

Lysozyme 4mg/ml, pI 11.3


Total Protein Loaded: 5.175mg

**Reproducible
separation over
target ionic capacity
range**

life
technologies

POROS[®] XS Characterization

Protein separation across the Ionic Capacity range

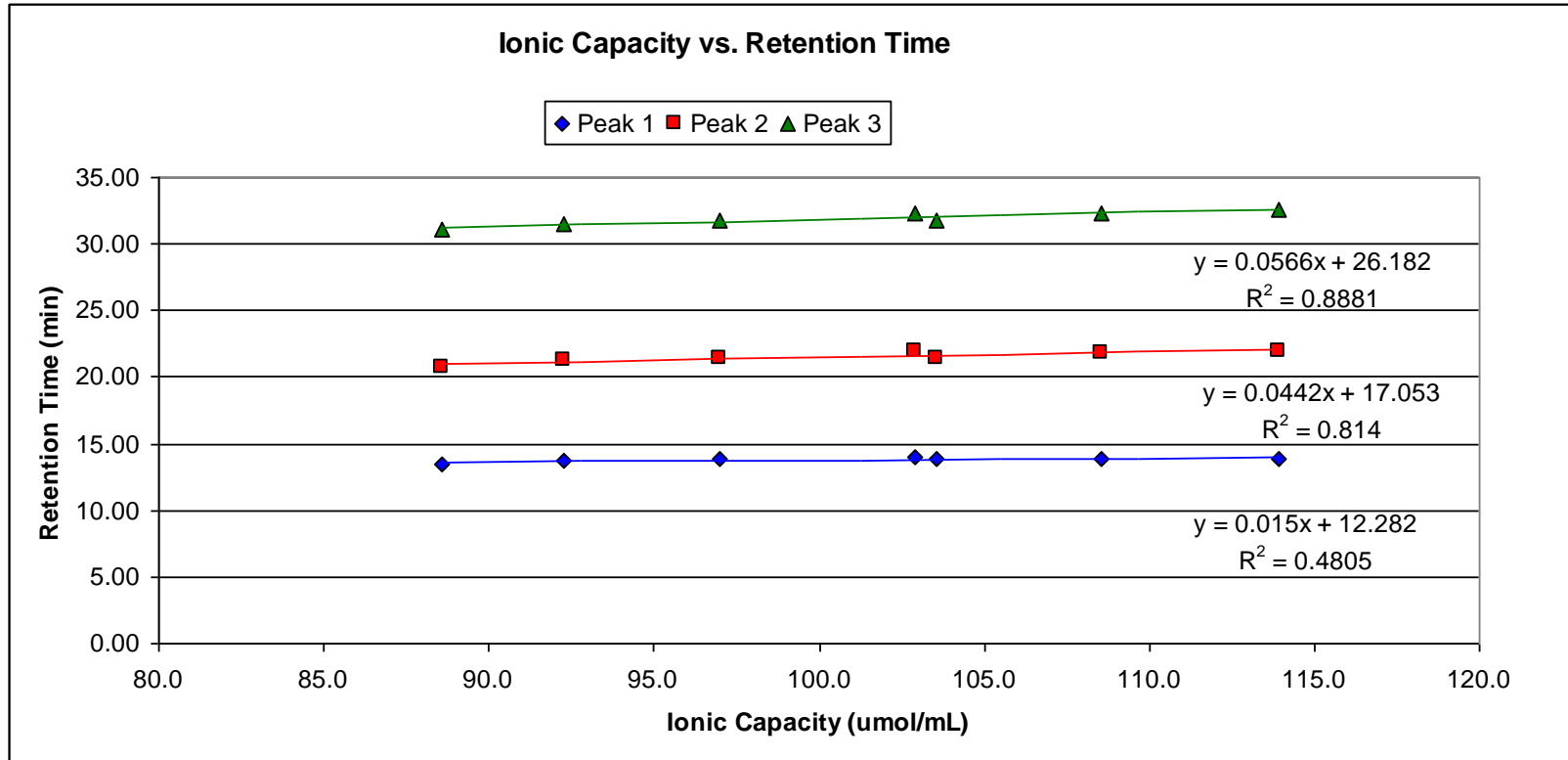


Minimal change in retention time
over target ionic capacity range



POROS[®] XS Characterization

Separation vs Protein Load

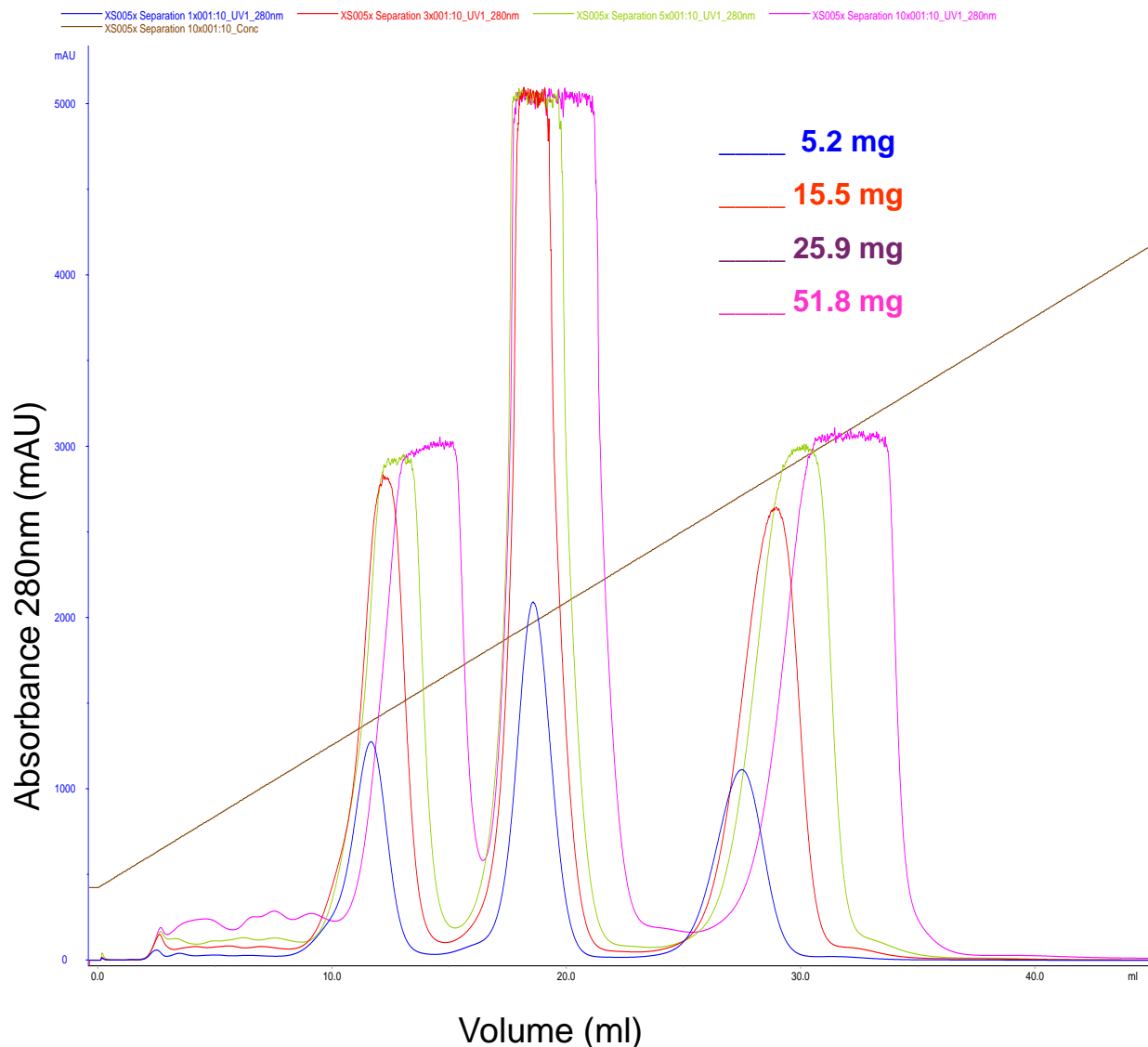

Format


Instrument: AKTA Purifier 10
Column: 0.46cmDx 20cmL, 3.32 ml
Gradient: 10 – 100%B, 16 CV
Buffer A: 20 mM MES, 25mM NaCl pH 6.2
Buffer B: 20 mM MES/ 1.0M NaCl, pH 6.2
Sample:
Chymotrypsinogen A 5mg/mL, pI 9.4
Cytochrome C 8mg/mL, pI 10.4
Lysozyme 4.25mg/ml, pI 11.3

Excellent resolution maintained even as load concentration is increased 10X



POROS[®] XS Characterization

Separation vs Bed Height/Residence Time

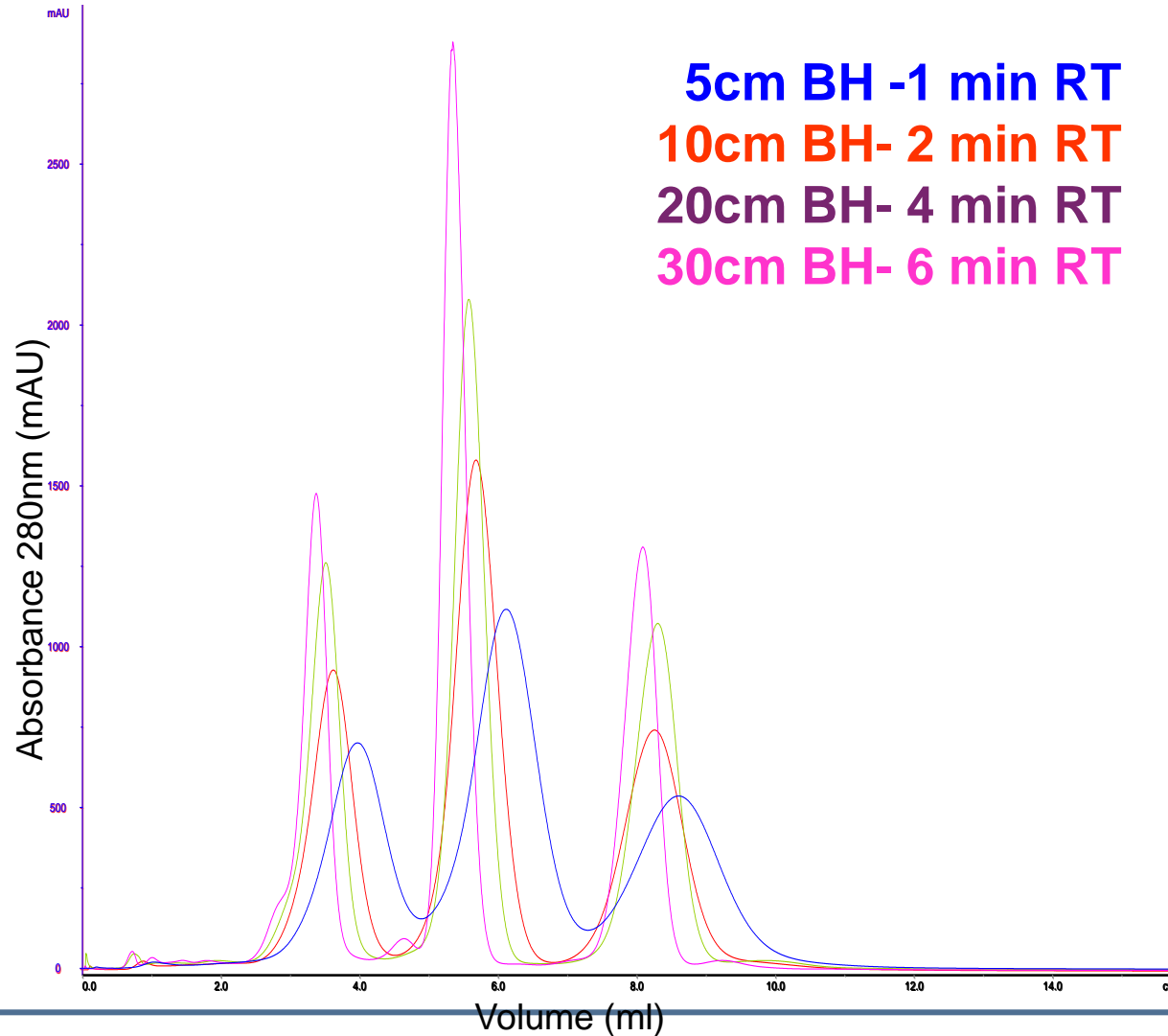

Format

Instrument: AKTA Purifier 10

Column: 0.46cmDx 20cmL, 3.32 ml

Gradient: 10 – 100%B, 16 CV

Buffer A: 20 mM MES, 25mM NaCl pH 6.2

Buffer B: 20 mM MES/ 1.0M NaCl, pH 6.2

Sample:

Chymotrypsinogen A 5mg/mL, pI 9.4

Cytochrome C 8mg/mL, pI 10.4

Lysozyme 4.25mg/ml, pI 11.3


**Small scale results
predict larger scale
allowing for flexible
scalability during
design**

POROS® XS Chromatography Resin

Setting Standards in Resolution and Capacity


Capacity of Cation Exchange Resins

40 mM NaCl, pH 5.0, 300 cm/H

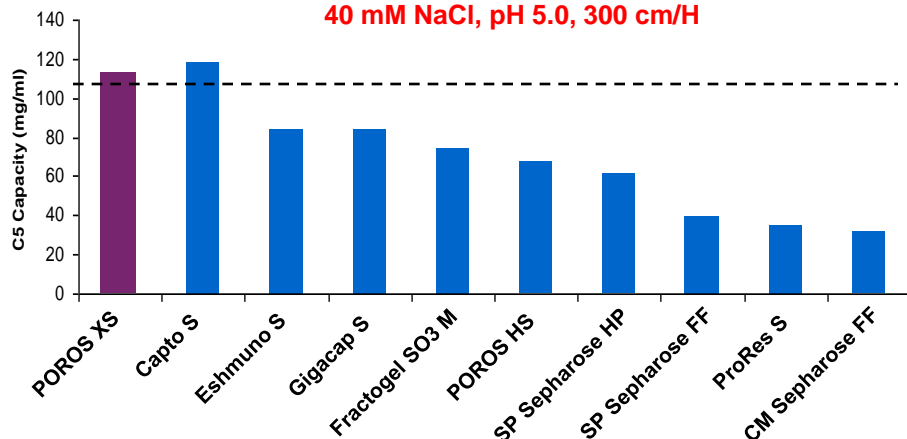

Capacity of Cation Exchange Resins

100 mM NaCl, pH 4.5, 300 cm/H

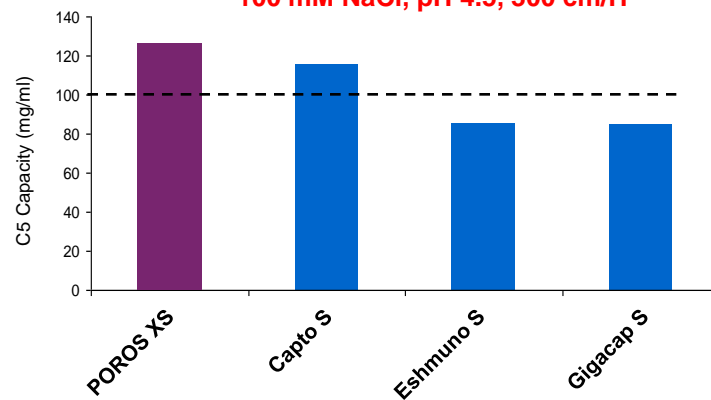


POROS® XS is the only CEX resin to successfully combine superior capacity, salt tolerance & resolution


Format

Column: 1cmD x 20cmL; **Buffer A:** 20mM MES, 25mM NaCl pH 6.2; **Buffer B:** 20mM MES, 1M NaCl pH 6.2;
Gradient: 10% B – 50% B 7.5CV; **Flow Rate:** 300cm/H; **Sample:** Chymotrypsinogen; Cytochrome C; Lysozyme


POROS XS vs. GigaCap S

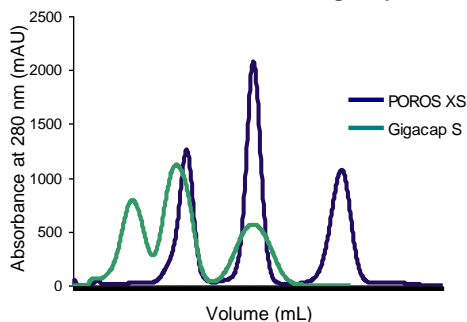

POROS XS vs. SP Sepharose FF

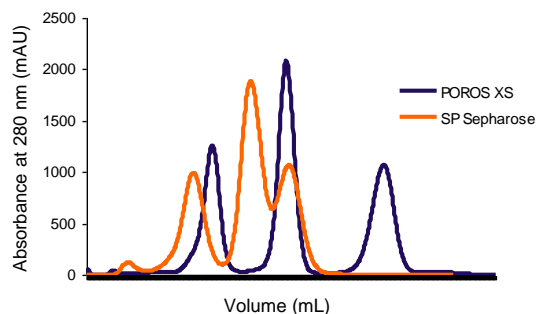

POROS XS vs. Capto S

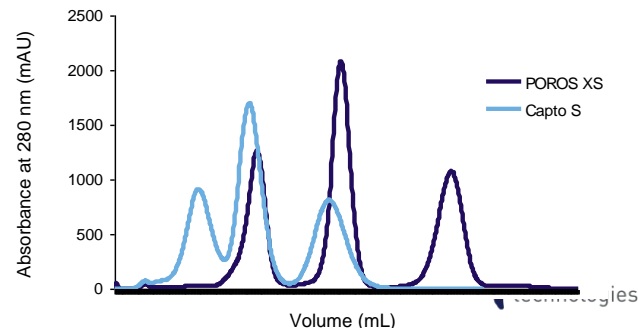


POROS[®] XS Chromatography Resin

Resolution


Superior resolution capabilities include:

- Maintaining the superior resolution of POROS[®] HS
- Providing consistent resolution across a wide range of flow rates
- Retaining excellent resolution through a ten-fold increase in load


POROS[®] XS and POROS[®] HS deliver equivalent resolution performance

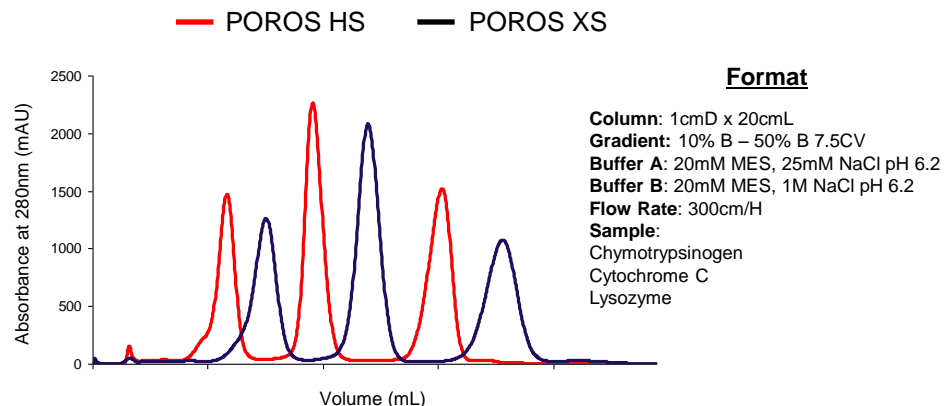


Excellent resolution is maintained as linear flow rate is increased.

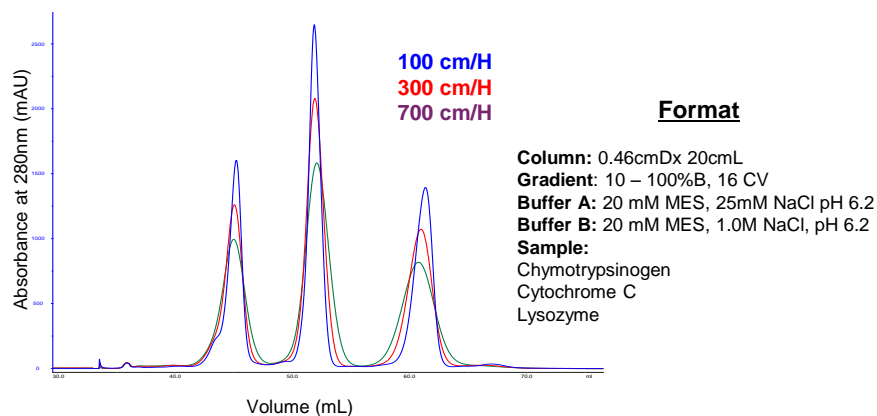

Resolution is maintained even as load concentration is increased 10x.

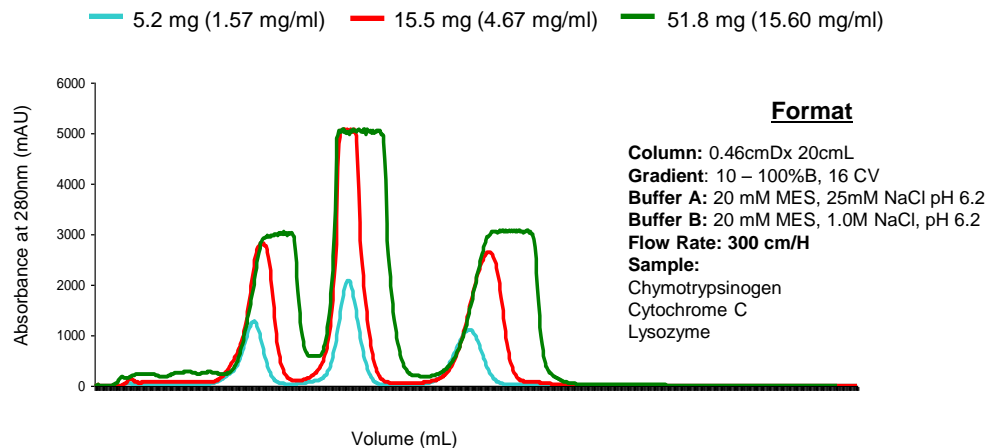

POROS[®] XS IgG binding capacity vs residence time

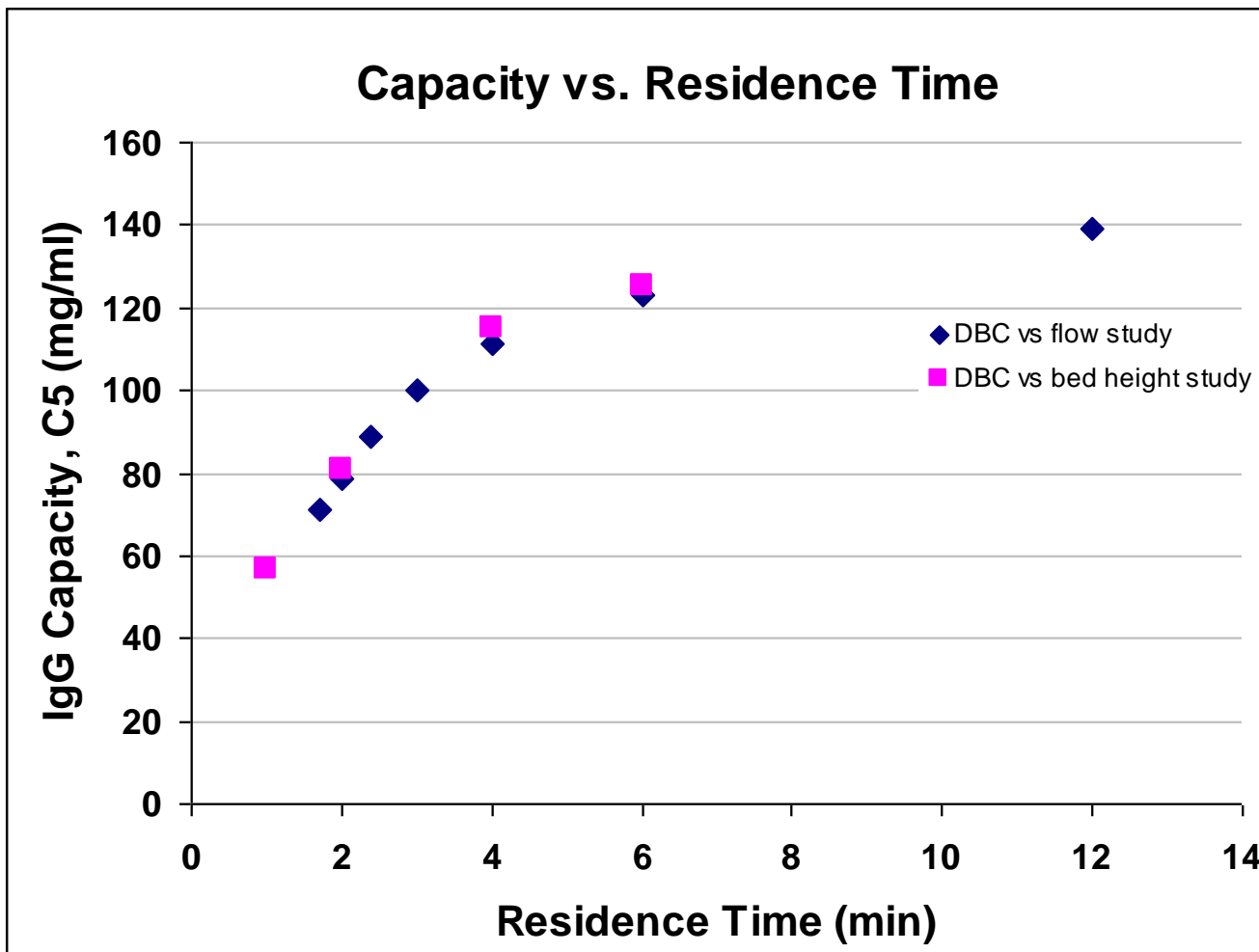


Format
Column Diameter: 0.46cmD
Buffer: 20mM MES, 40mM NaCl pH 5.0
Load Condition: 5 mg/ml Polyclonal IgG in Equil

Binding capacity is ≥ 100 mg/ml at a residence time of 3 minutes and greater



POROS[®] XS IgG binding capacity vs bed height

POROS XS: IgG Capacity vs Bed Height

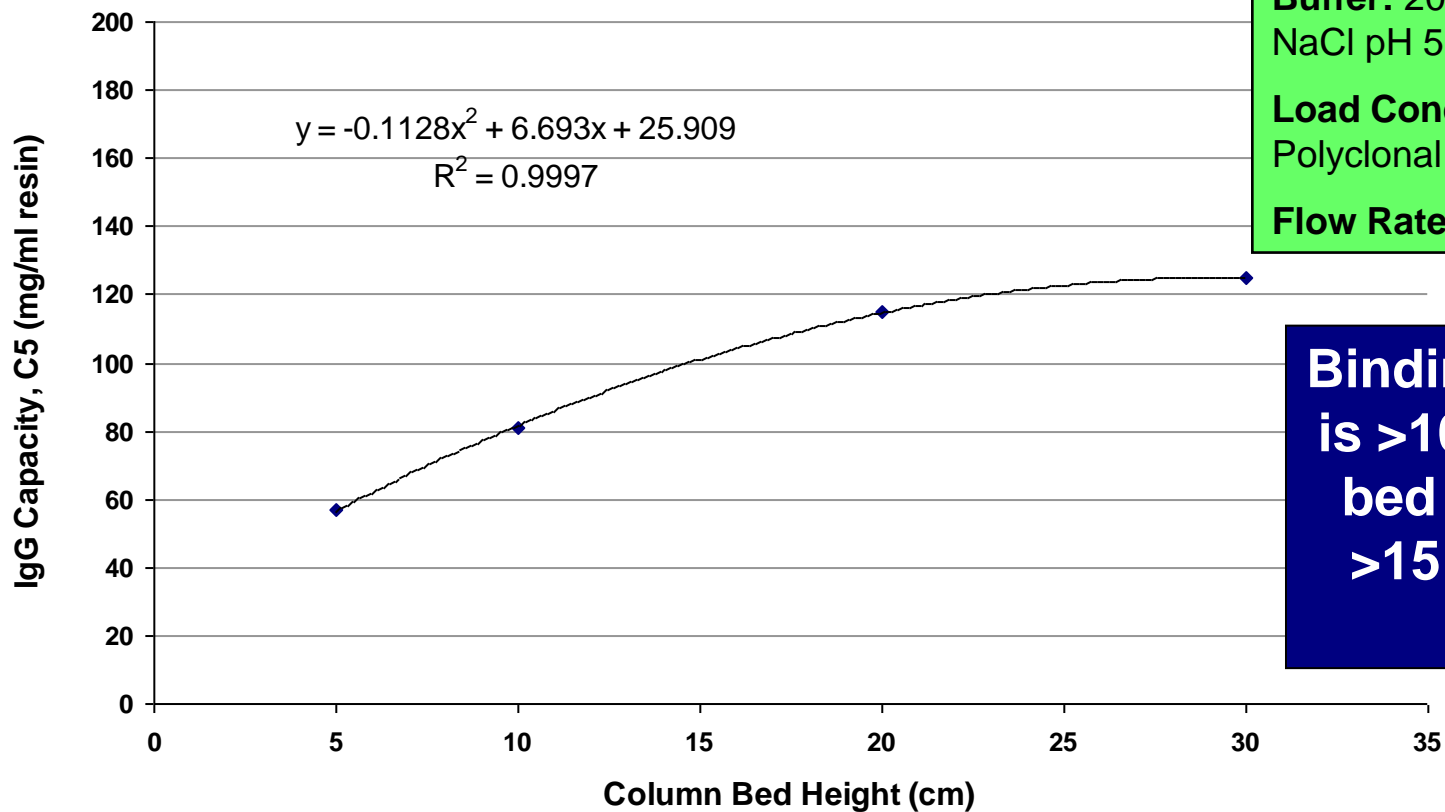

Format

Column Diameter: 0.46cmD

Buffer: 20mM MES, 40mM NaCl pH 5.0


Load Condition: 5 mg/ml Polyclonal IgG in Equil

Flow Rate: 300 cm/hr

Binding Capacity is >100 mg/ml at bed heights of >15 cm at 300 cm/hr



POROS[®] XS IgG capacity vs load concentration

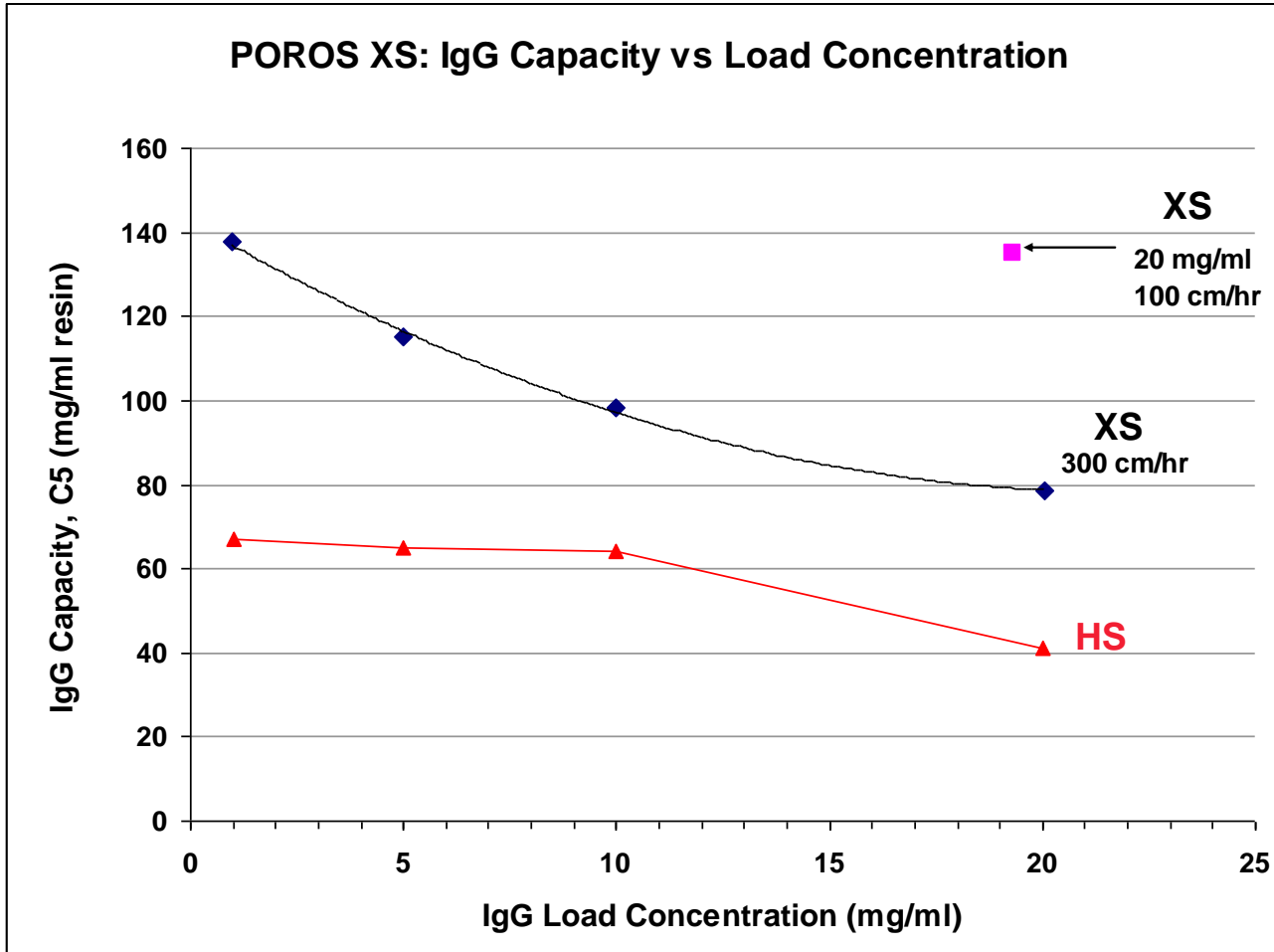


Format
Column: 0.46cmDx 20cmL
Buffer: 20mM MES, 40mM NaCl pH 5.0
Load Condition: Polyclonal IgG in Equil
Flow Rate: 300 cm/hr

Binding Capacity is >80 mg/ml even at concentrations >20 mg/ml



POROS[®] XS IgG Binding Capacity vs Buffer System

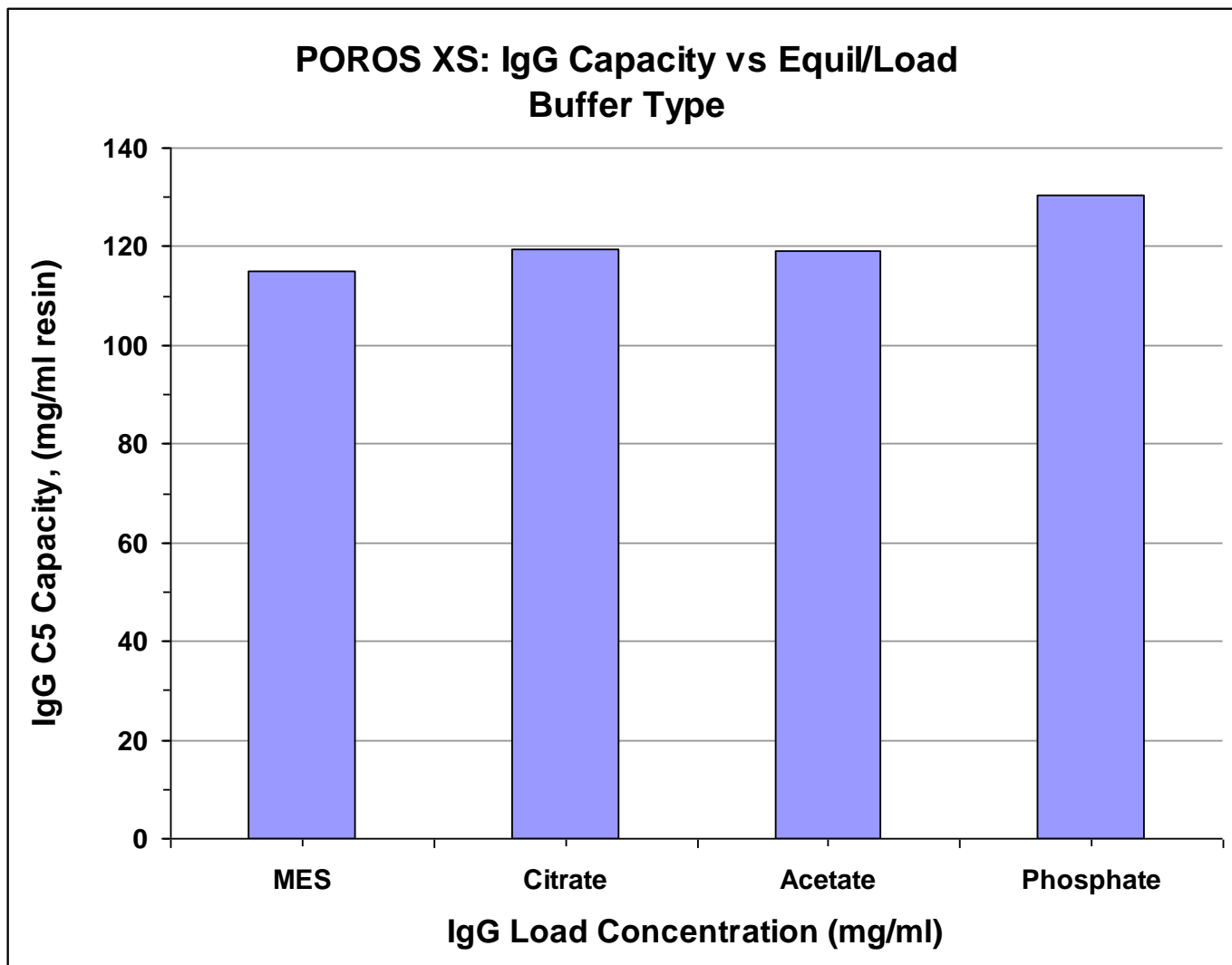


Format

Column: 0.46cmDx 20cmL
Buffer: 20mM MES, citrate, acetate, phosphate
Cond: 4.7ms/cm
pH: 5.0
Load Condition: 5 mg/ml Polyclonal IgG in Equil
Flow Rate: 300 cm/hr

The buffer system used does not impact binding capacity



POROS[®] XS IgG Binding Capacity vs Salt Concentration

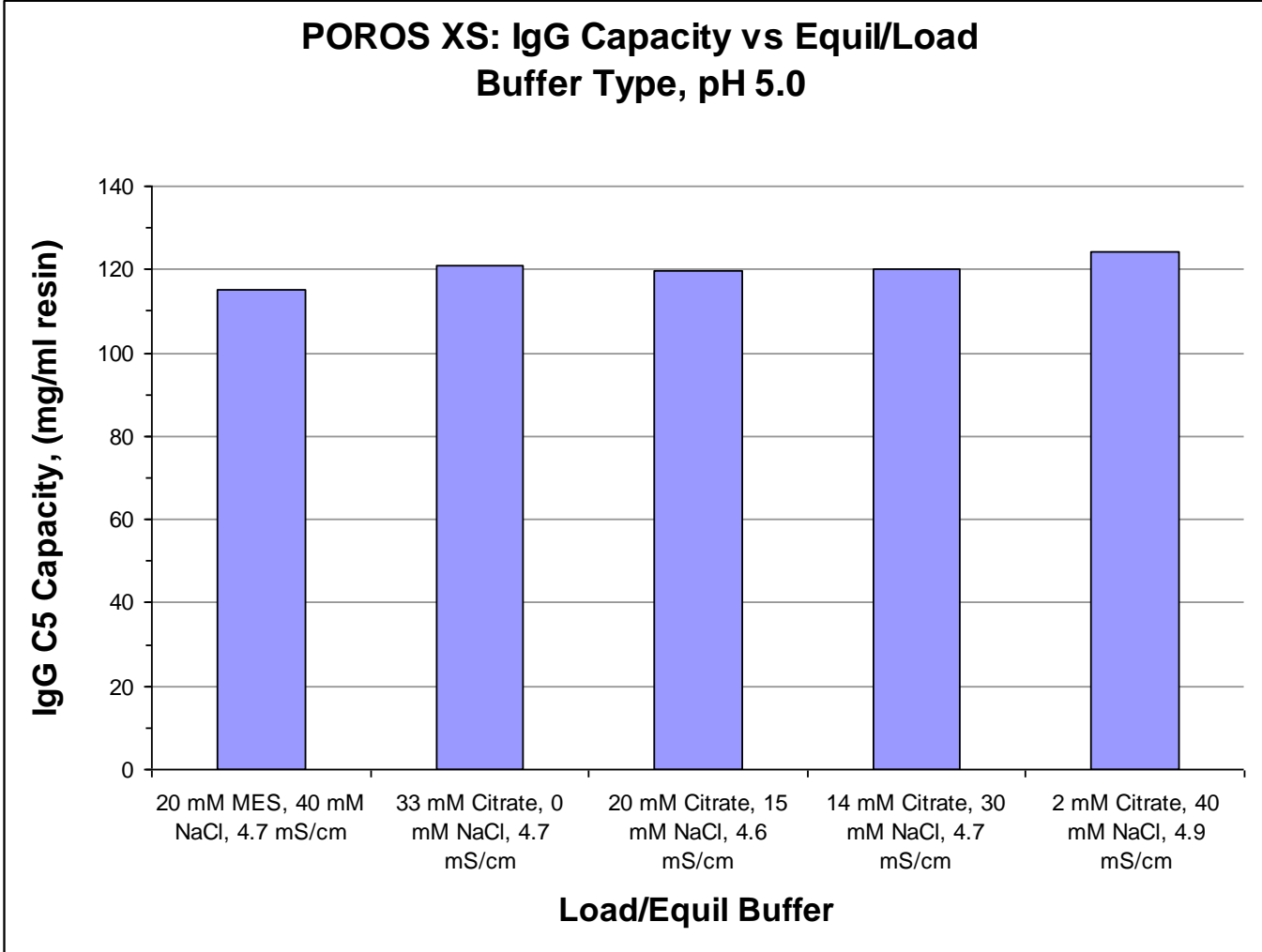

Format
Column: 0.46cmDx 20cmL
Buffer: citrate/NaCl in varying ratios
Cond: 4.7ms/cm
pH: 5.0
Load Condition: 5 mg/ml Polyclonal IgG in Equil
Flow Rate: 300 cm/hr

Binding capacity not affected by NaCl concentration at a given conductivity



POROS® XS Protein Separation vs. Buffer System

Separation Control Conditions

Column = 0.46 cmD x 20 cmL


Buffer A = 20mM MES, 25mM NaCl
pH 6.2 C = 4.02 mS/cm

Buffer B = 20mM MES, 1M NaCl pH
6.2 C = 89.7 mS/cm

Elution = 0% B – 100% B 16.88 CV

Flow Rate = 300cm/H

Protein Mix = Chymotrypsinogen A
Cytochrome C
Lysozyme




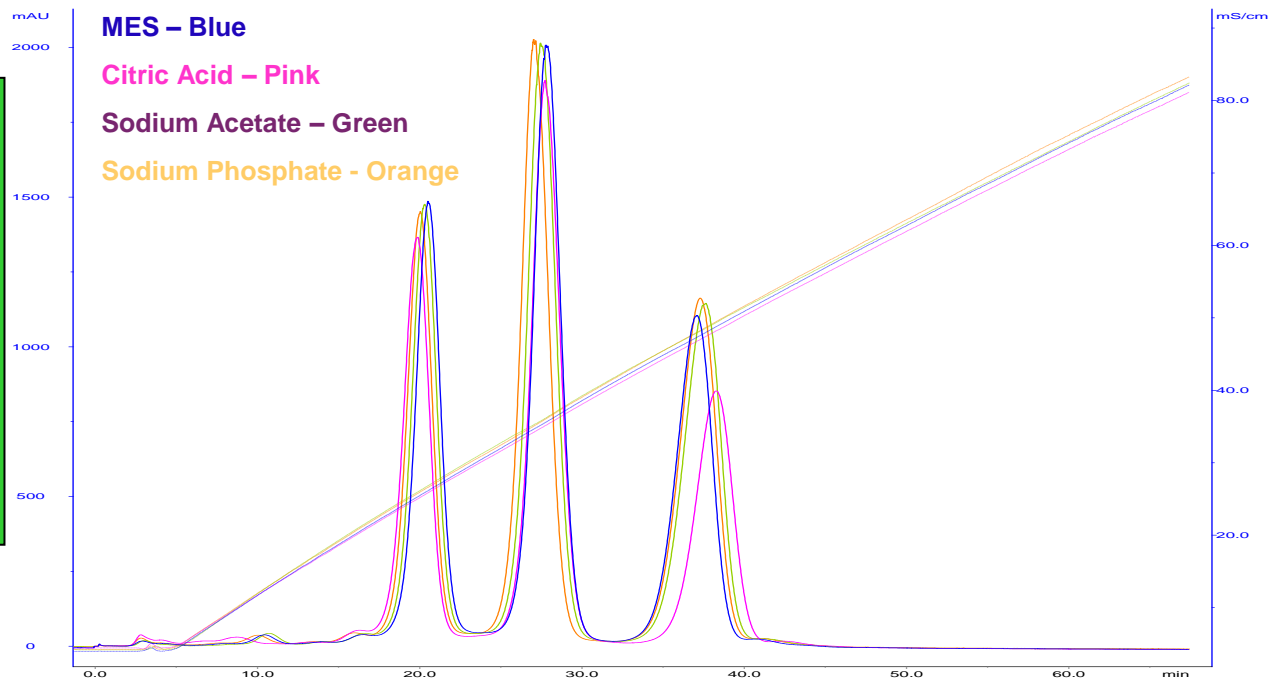


Sample	Retention Time (min)			Elution Conductivity (mS/cm)			Resolution	
	Peak 1	Peak 2	Peak 3	Peak 1	Peak 2	Peak 3	Peak 1 - 2	Peak 2 - 3
MES	20.50	27.77	37.10	26.120	35.589	47.299	2.45	2.55
Citric Acid	19.89	27.73	38.28	25.004	35.112	48.178	2.52	2.80
Sodium Acetate	20.29	27.43	37.63	26.517	35.742	48.568	2.36	2.68
Sodium Phosphate	20.03	27.01	37.27	25.981	35.044	48.261	2.33	2.76

POROS® XS maintains excellent resolution with multiple buffer conditions

Overview

- CEX introduction: POROS® HS and XS
- Applications and Characterization Data
 - Applications
 - Viral clearance
 - Pressure vs flow curves
 - Particle size analysis
- How to optimize
- Case Studies
 - POROS® HS process redesign
 - POROS® HS in HTP screening and scale up study
 - Rethinking polish applications: CEX in flow through/overload mode
- Summary

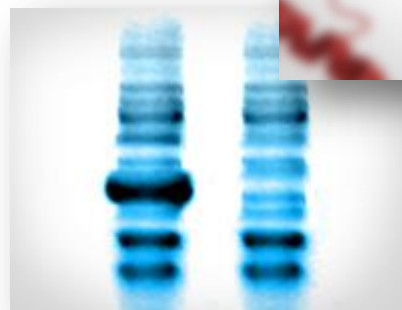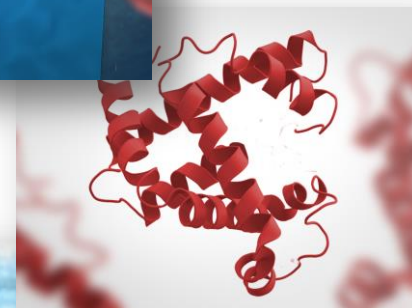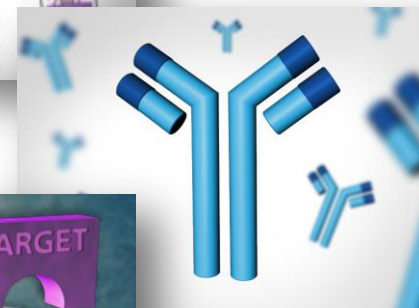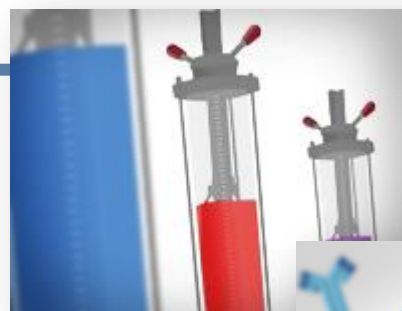


Optimization of a Conventional Soft Gel CEX Resin

Run 1 – pH 5.5

Run 2 – pH 5.0

Run 3 – pH 4.5

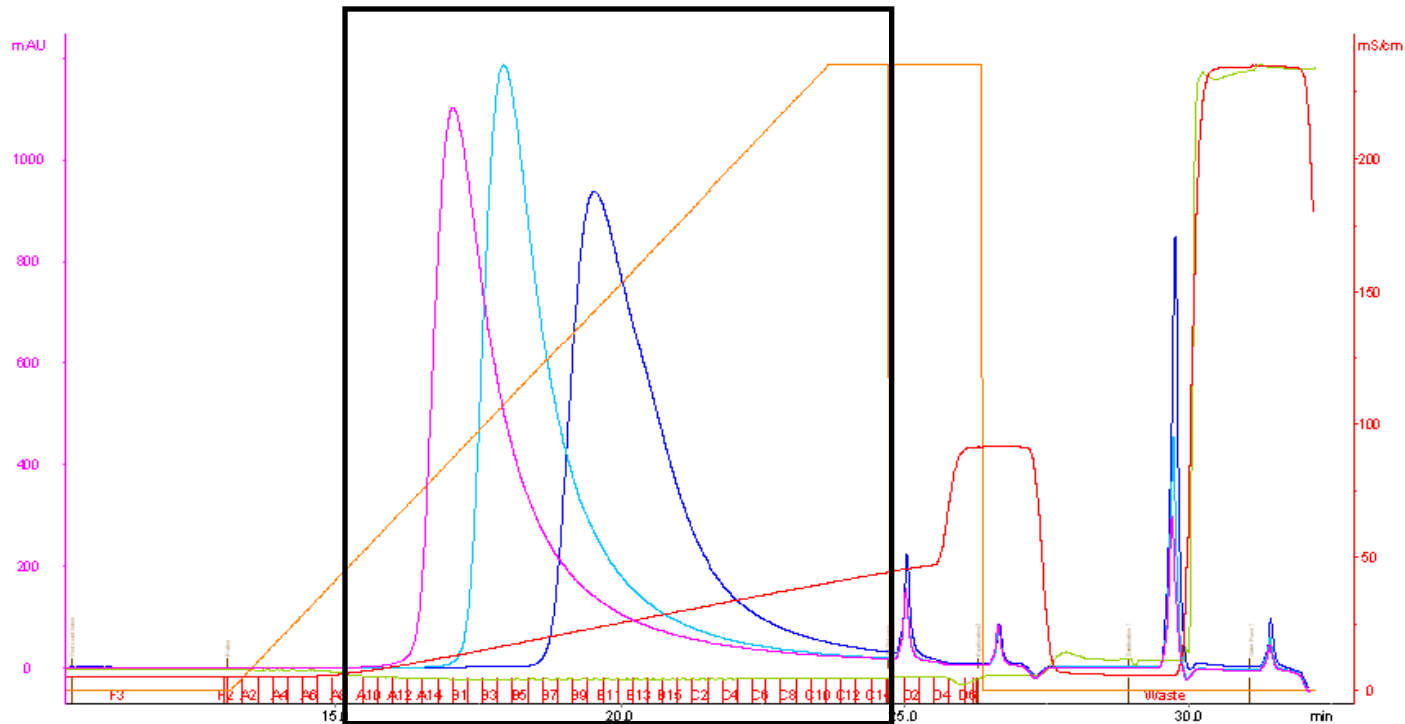


- Large particle size resin lacks resolution capability
- Peak cutting required to separate monomer and aggregate
- Decreased product yield

Optimization on POROS[®] XS

Run 1 – pH 5.5


Run 2 – pH 5.0

Run 3 – pH 4.5

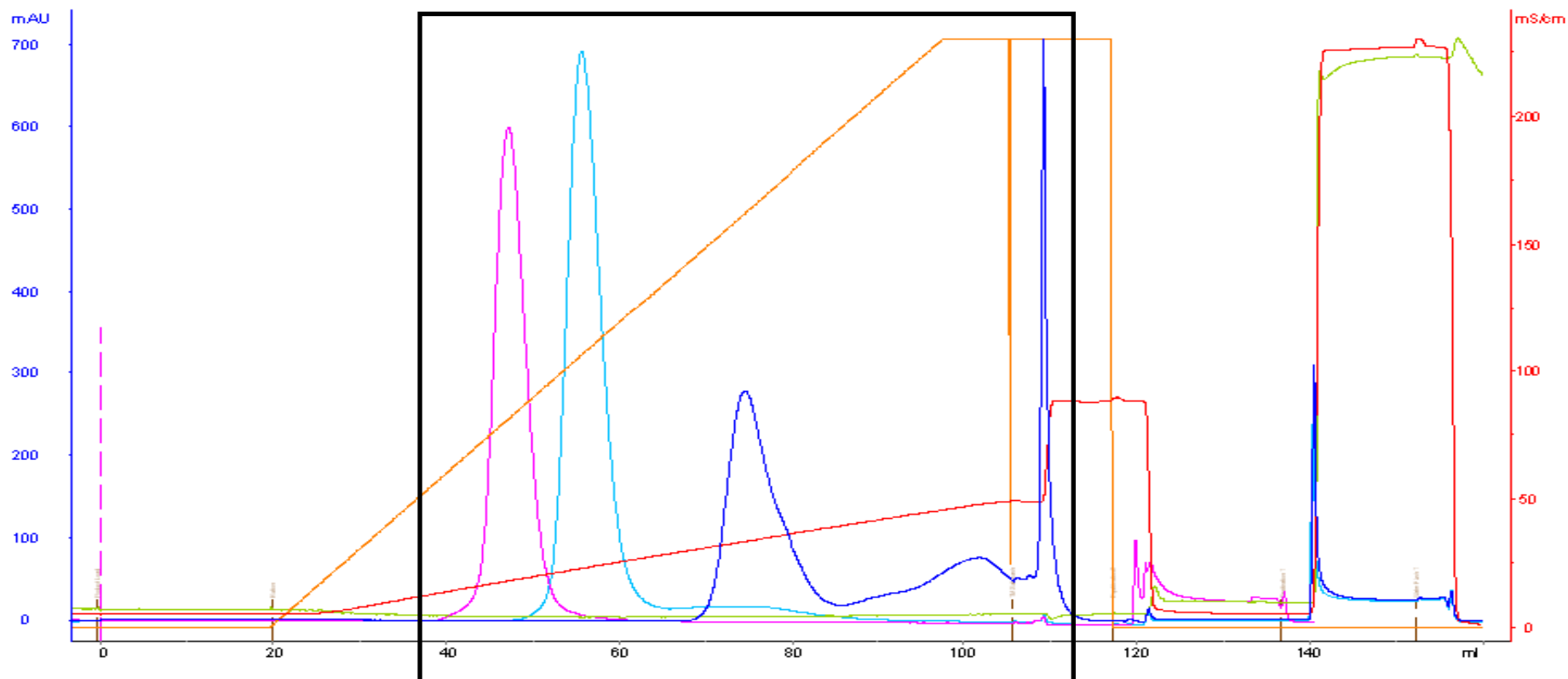

High resolution POROS[®] Resin provides improved separation

Optimization on POROS® XS

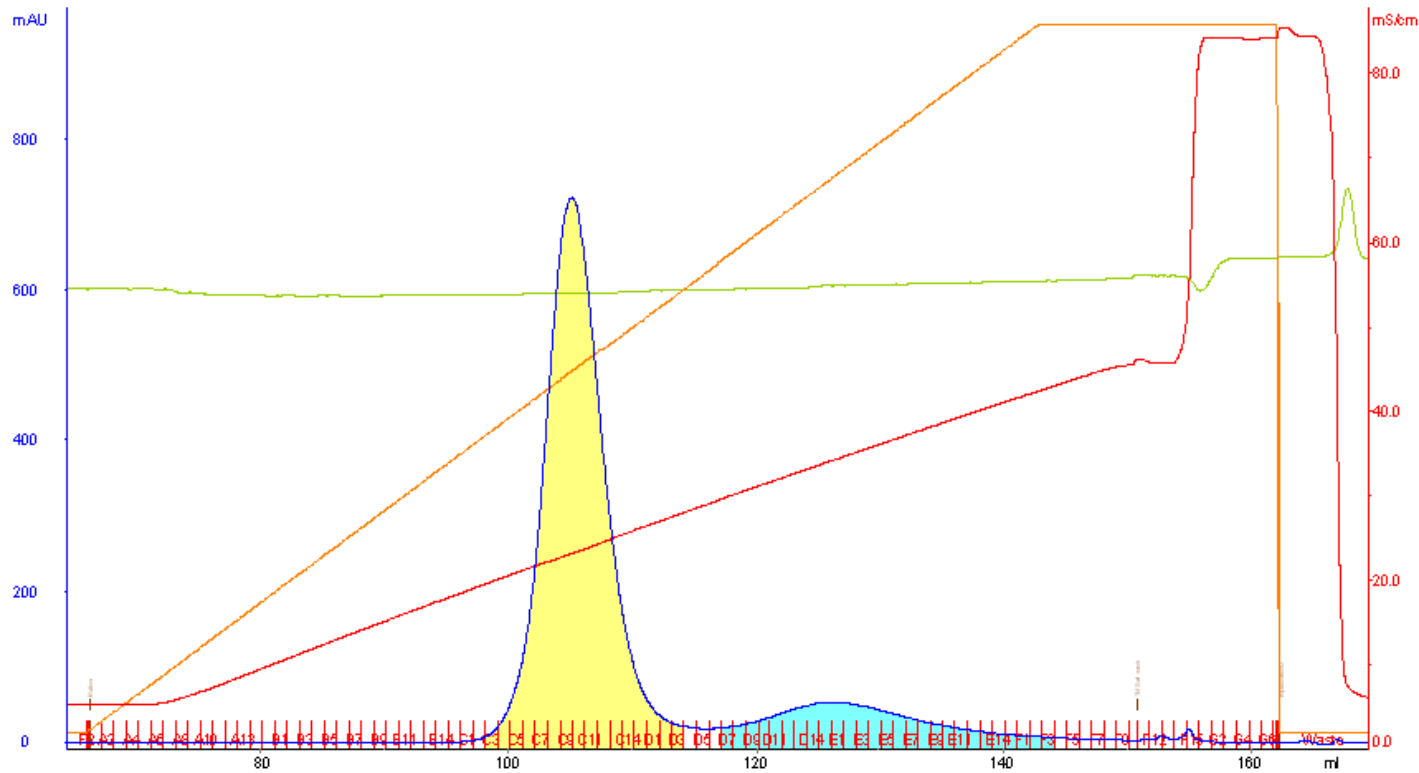

- Optimization of baseline separation
- Gradient can be optimized to step elution
- Improved product yield due to efficient separation of monomer and aggregate

POROS[®] XS Elution Peak: Customer Application

Goal: Ability to remove leached protein A below 10 ppm and maintain high process yield and product purity

Format

Column = 0.46cmD x 20cmL


Load : 70 mg/ml of resin, ≤ 5mS/cm, pH 5.0 by adjusting the neutralized Protein A elution peak with 1M Citric Acid and water

Gradient: 0 – 100%B, 20 CV

Buffer A: 20mM Sodium Citrate pH 5.0

Buffer B: 20mM Sodium Citrate, 500mM Sodium Chloride pH 5.0

Fractions: ½ CV, 1.662mL

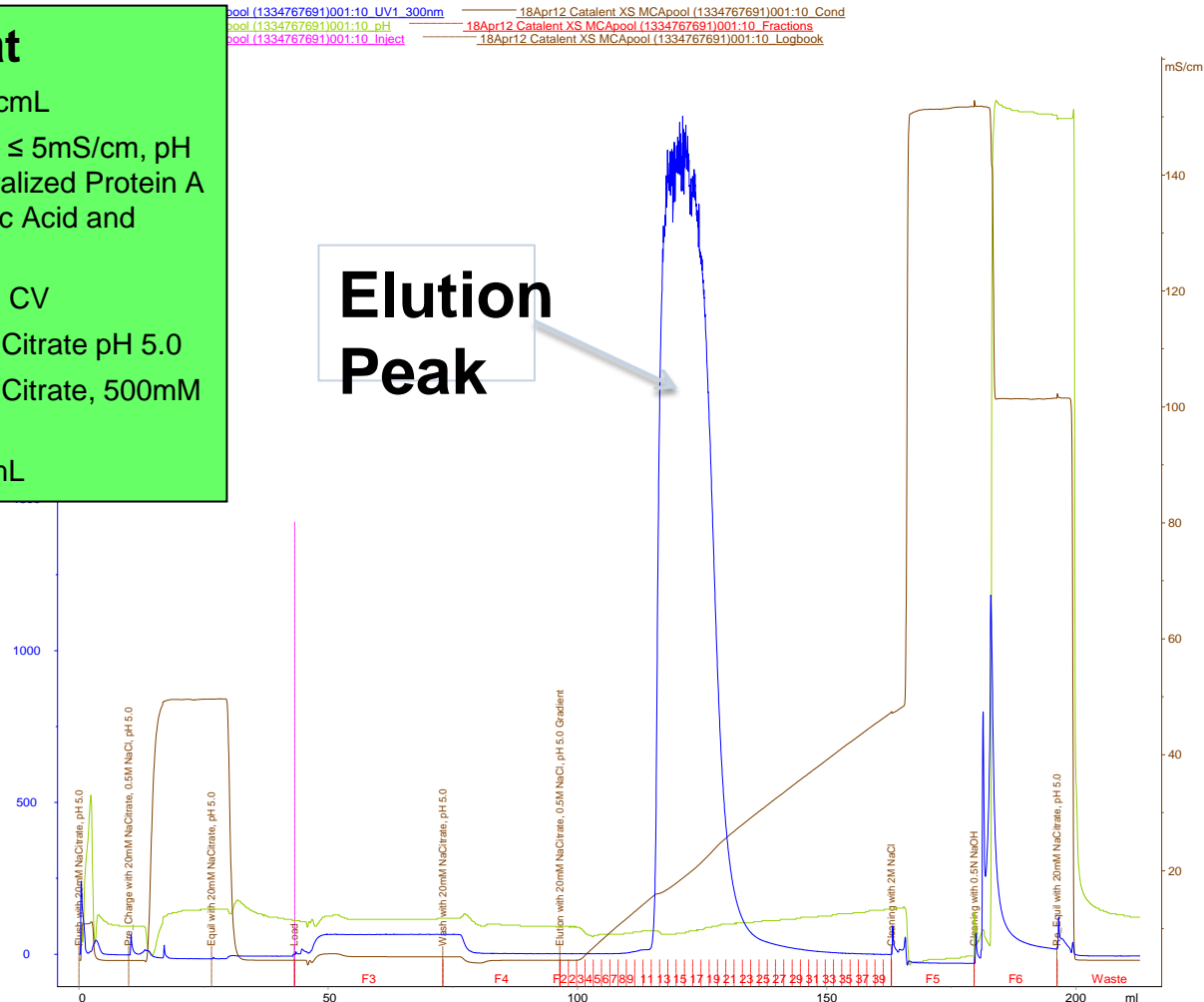


life
technologies

POROS® XS Product Pool

Yield & Purity

Sample	% HMW	Protein A (ppm)	% Yield
XS Load	1.007	40.65	
XS Elution Pool	0.533	Below LOD	94.5%

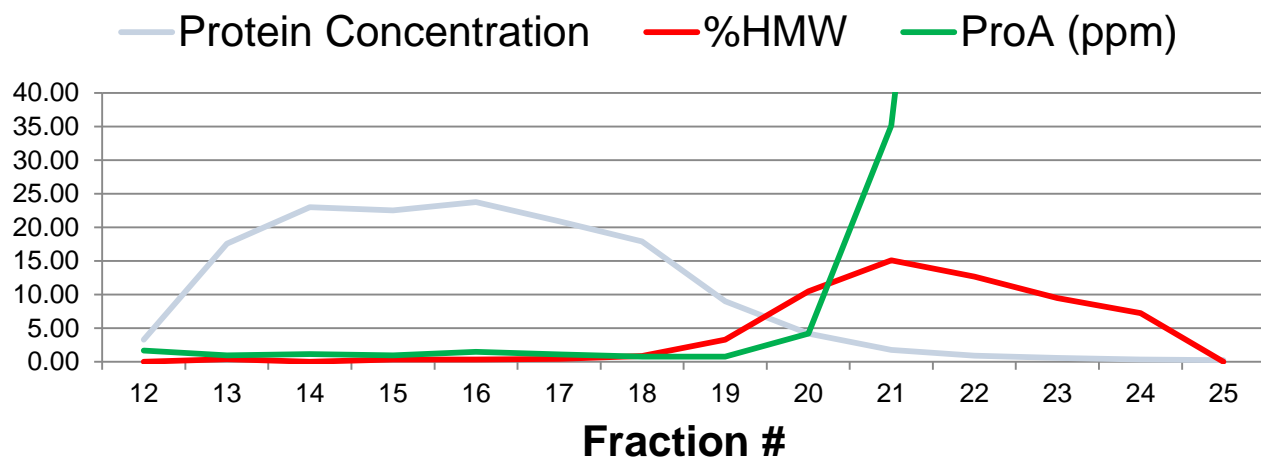

- High resolution yields effective separation of aggregate and leached rPA with high monomer yield
- SEC-HPLC analysis showed aggregate levels were reduced 2-fold
- Leached protein A reduced to below LOD
- Leached protein A was well separated from the MAbs peak (salt strip) and final CEX yield was high

Viral Clearance on POROS XS in Bind/Elute Mode

XMuLV Clearance (Log10)

Load pH	pH 6.0	pH 5.0	pH 4.5
Load Salt Concentration (mM)	25	25	100
Load Capacity (mg/ml)	80	80	100
FT/Wash	3.5	4.0	2.1
20 mM MES, 50 mM NaCl		3.0	
20 mM MES, 100 mM NaCl	1.8	1.8	
20 mM MES, 200 mM NaCl	1.1	1.9	
20 mM MES, 300 mM NaCl	1.0	2.0	3.1
20 mM MES, 400 mM NaCl		1.9	
20 mM MES, 500 mM NaCl	1.0	1.5	1.7
2M NaCl		1.1	

Study Conditions

Column Format: 0.46cmD x 20 cmL, 3.3 ml

Flow Rate: 300 cm/hr


Residence Time: 4 min

Load: 5 mg/ml Human Polyclonal IgG

Virus Spike: 5% Vol


- Salt concentration in the load appears to increase virus flow through on POROS XS (lower clearance in FT/Wash) more than pH in the range tested
- pH 4.5 with 100mM NaCl in load appears to have best clearance when eluting with 300 mM NaCl: 3.1 LRV
- POROS XS is capable of good viral clearance depending on the process conditions as with any cation exchange resin

POROS[®] 50 μ m: Pressure vs. Flow Curves

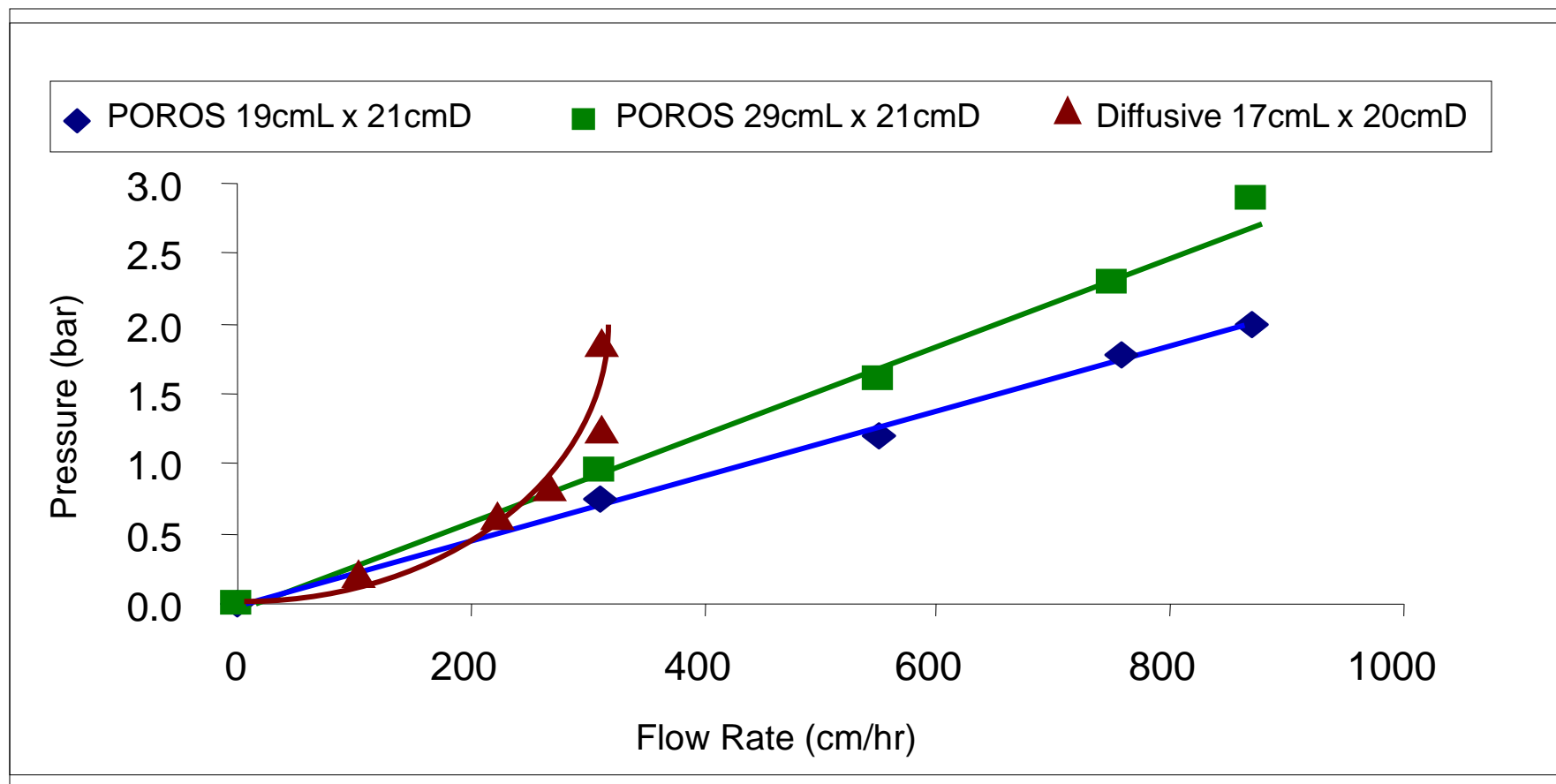

POROS[®] Resin: Linear pressure response with no bed compression


POROS[®] HS Chromatography Resin Pressure Flow 200 cm Diameter Column

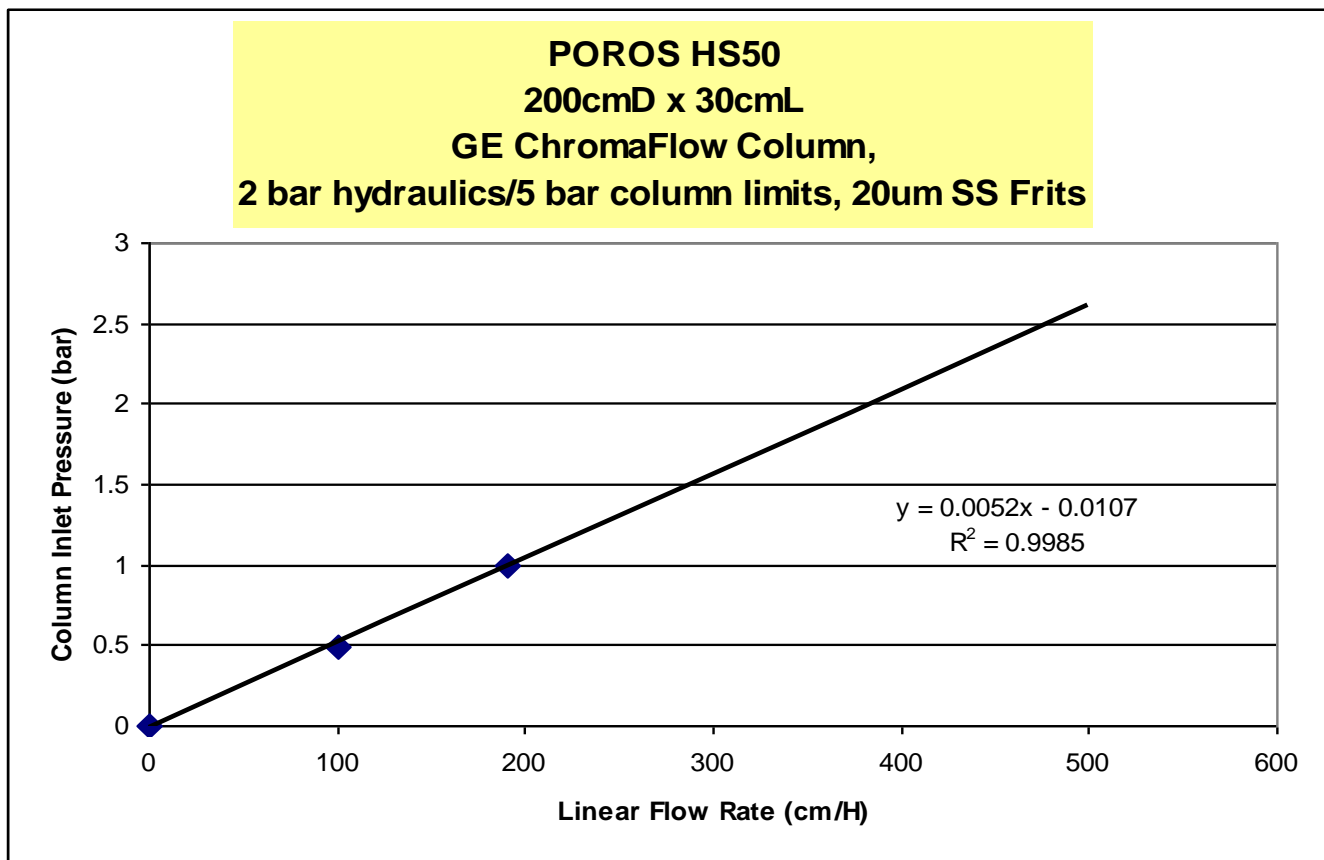

Note: The system pressure was not subtracted from data

POROS[®] XS Pressure vs Flow


6.2 cm Diameter Column

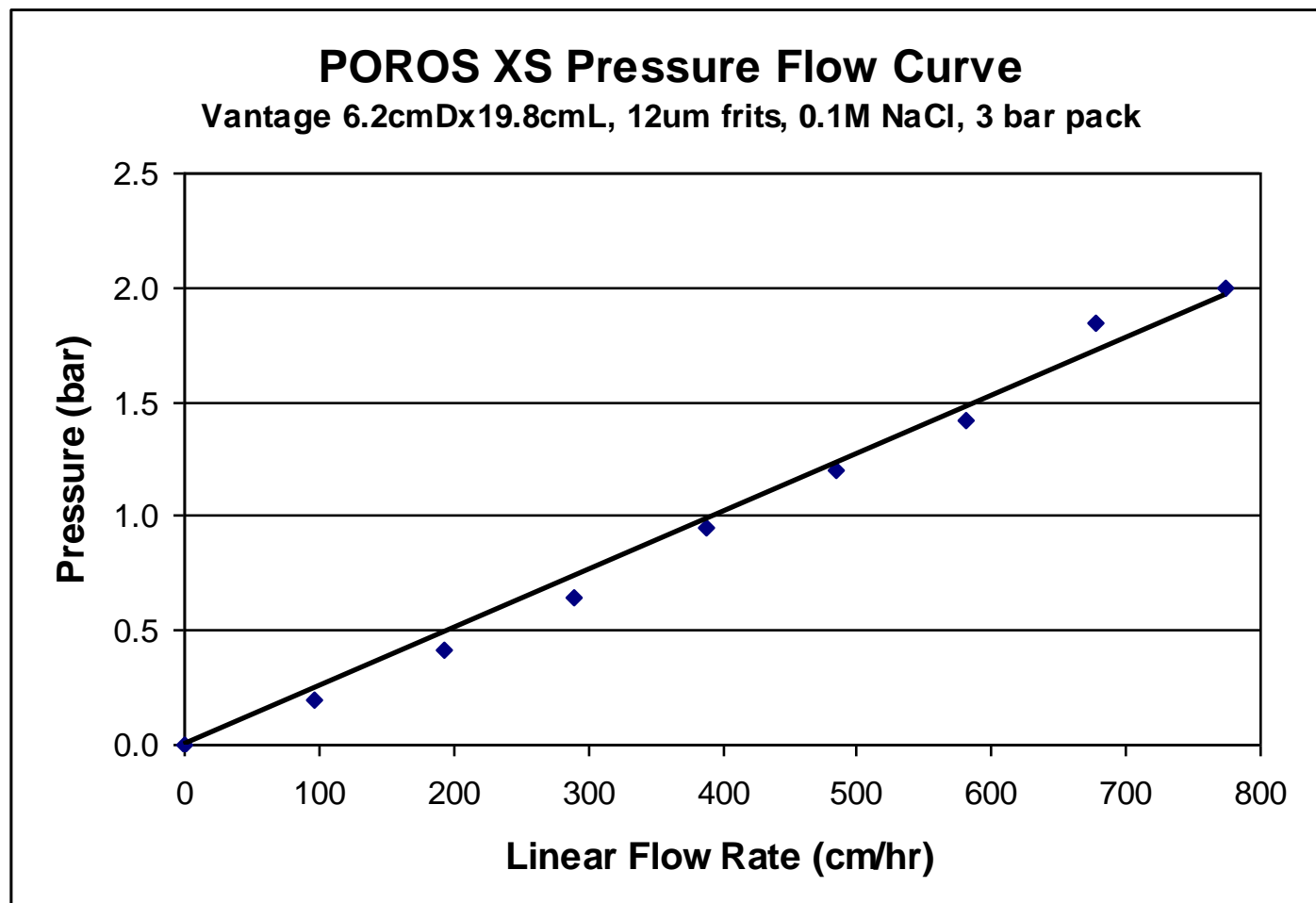

Linear pressure response

POROS[®] XS Pressure vs Flow

30 cm Diameter Column

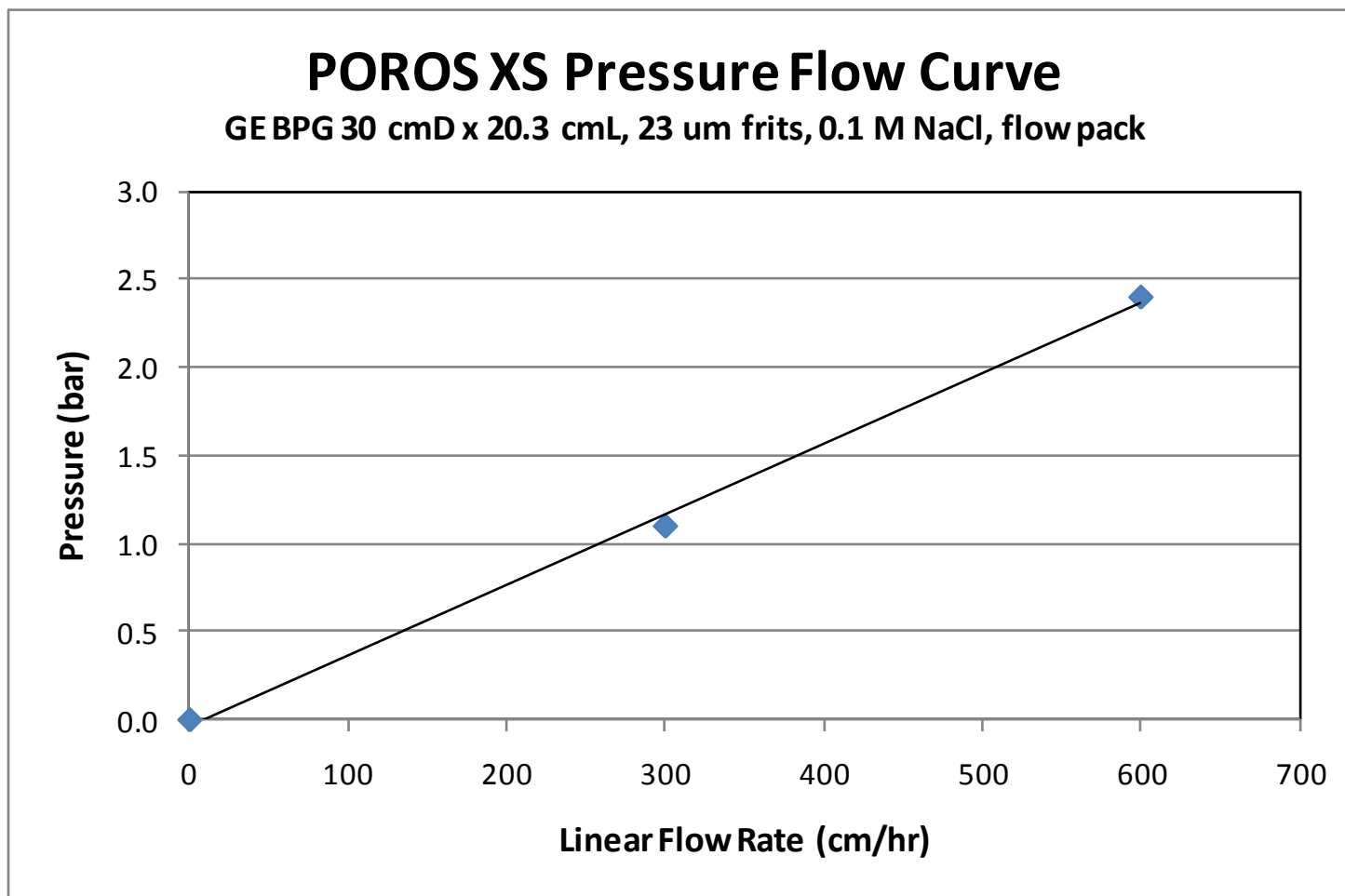


* System pressure not subtracted

life
technologies

Linear Pressure Response

POROS[®] XS and HS Pressure vs Flow

1 cm Diameter Column

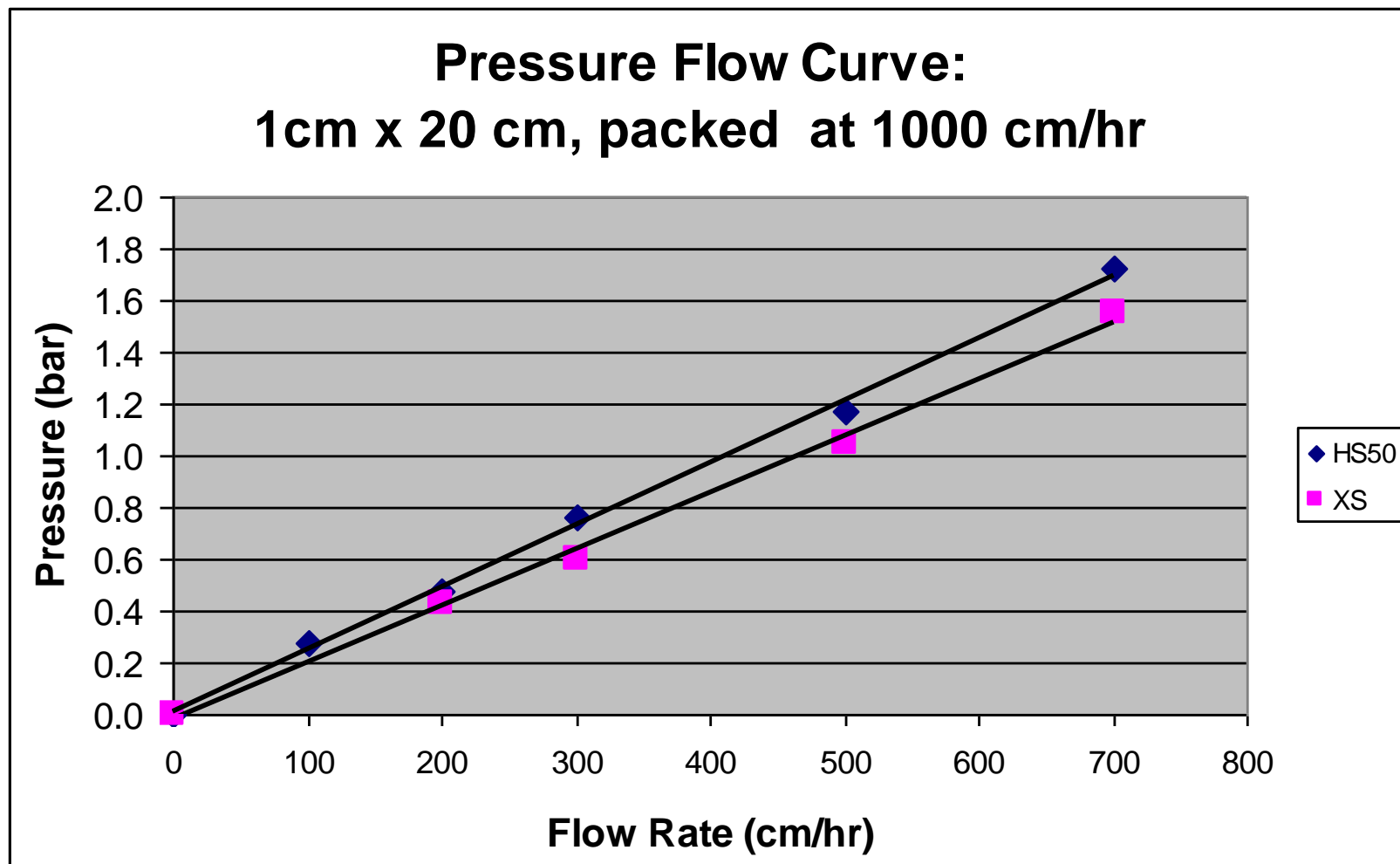


Linear pressure response

life
technologies

POROS[®] XS Pressure vs Flow

2.5 cm Diameter Column in 0.1M NaCl and water

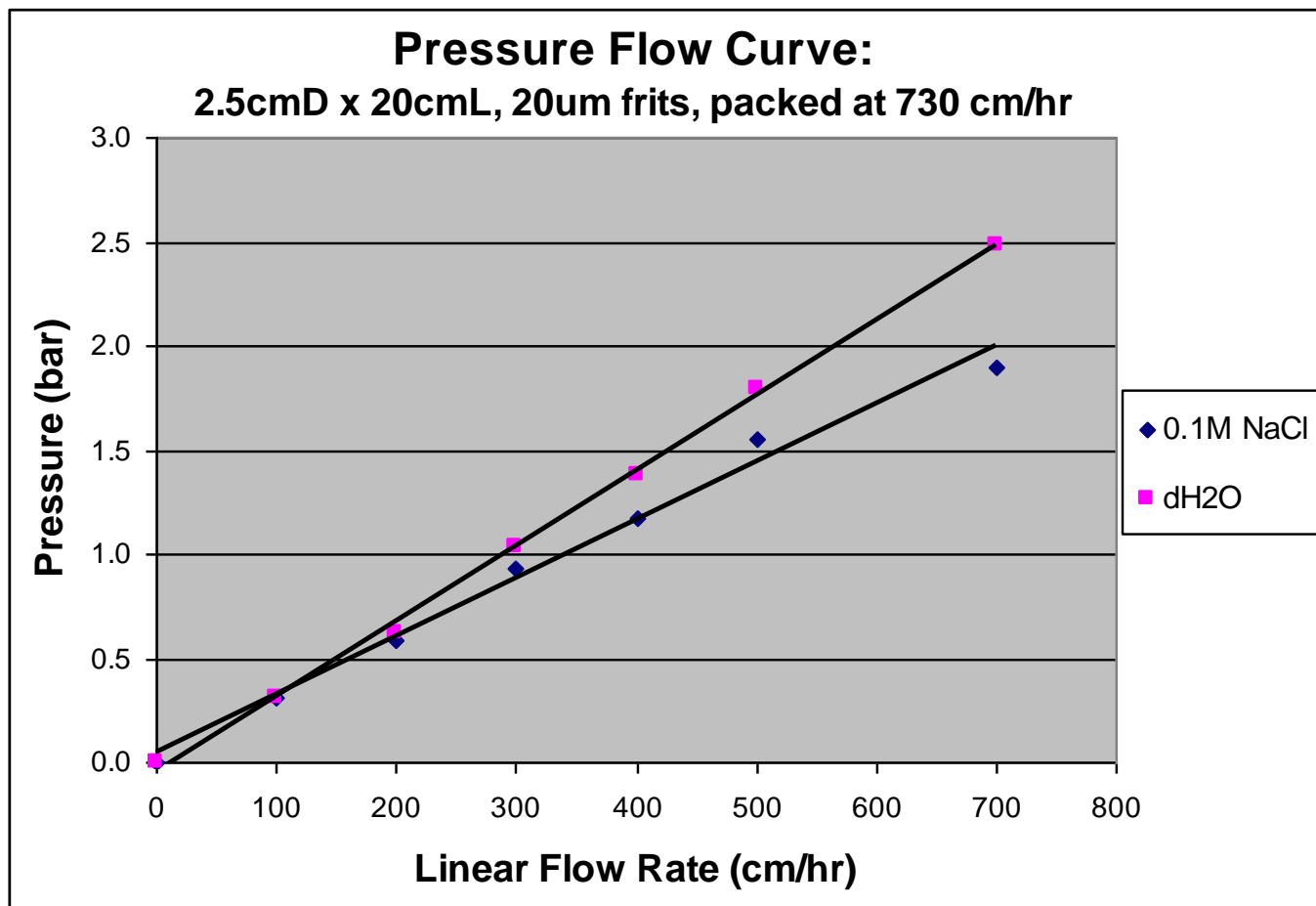


Linear pressure response

life
technologies

POROS[®] HS Chromatography Resin

Particle size distribution


- 50um Bulk Resin: Tight Particle Size Distribution

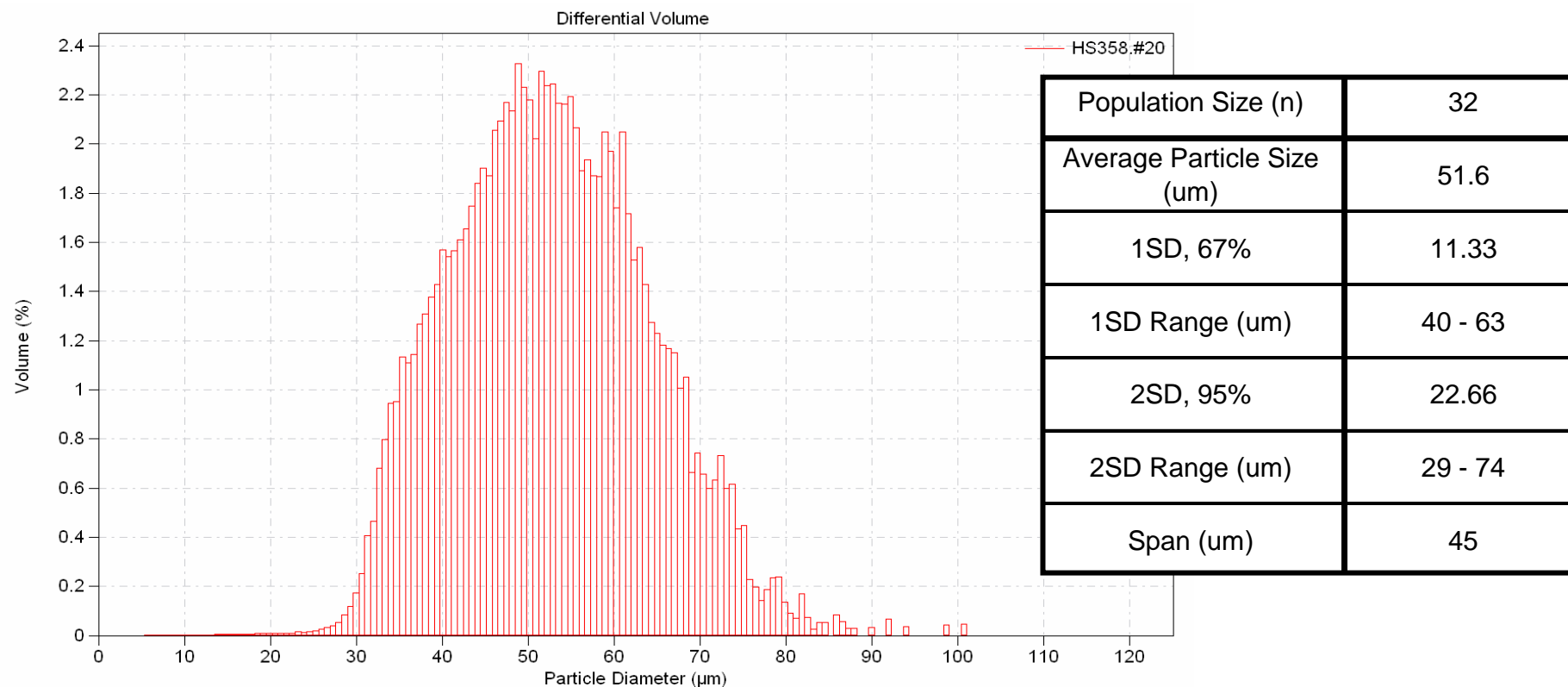

POROS[®] XS Chromatography Resin

Particle size distribution

- POROS[®] XS Resin: Tight Particle Size Distribution




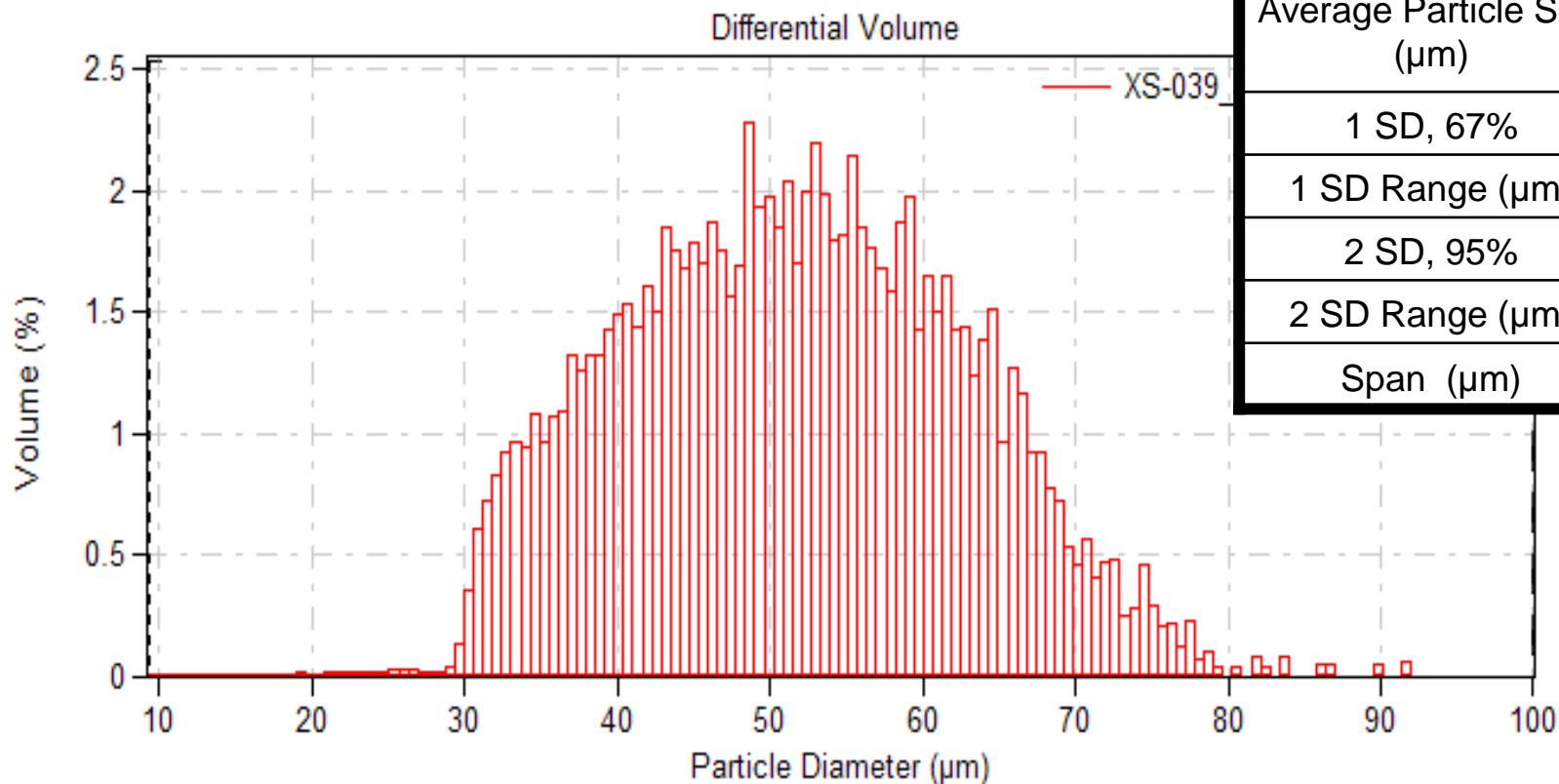


LC= 5.416 µm UC= 162.5 µm {100%}

life
technologies

Overview

- CEX introduction: POROS[®] HS and XS
- Applications and Characterization Data
 - Applications
 - Viral clearance
 - Pressure vs flow curves
 - Particle size analysis
- How to optimize
- Case Studies:
 - POROS[®] HS process redesign
 - POROS[®] HS in HTP screening and scale up study
 - Rethinking polish applications: CEX in flow through/overload mode
- Summary

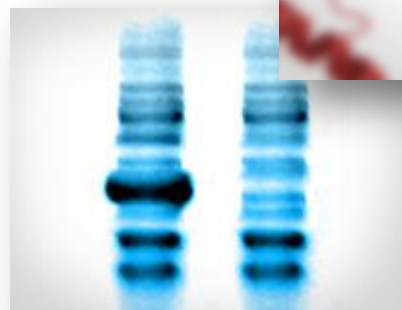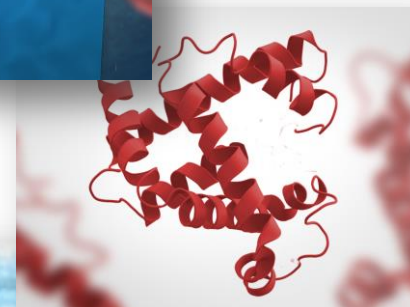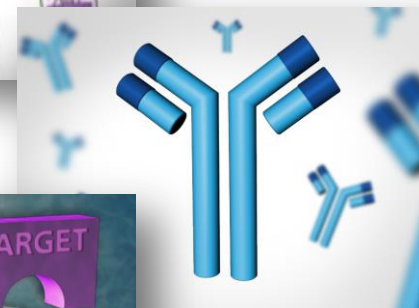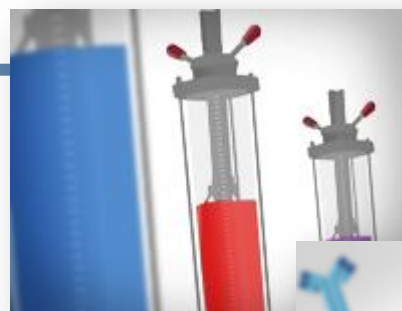


Optimizing POROS CEX Conditions: Where to start?

- Although similar sulphopropyl functional groups are used on most strong CEX resins, the optimal binding and elution conditions can vary significantly due to a number of resin characteristics
- Different CEX resins operated with the same process conditions will yield variable results; therefore, standardized conditions or platform-type evaluations are not recommended
- It is important to test different loading and elution conditions to optimize capacity, separation and yield based on the target molecule characteristics and process challenges

Benchmarking Cation Exchange Performance

IgG DBC surface response plots

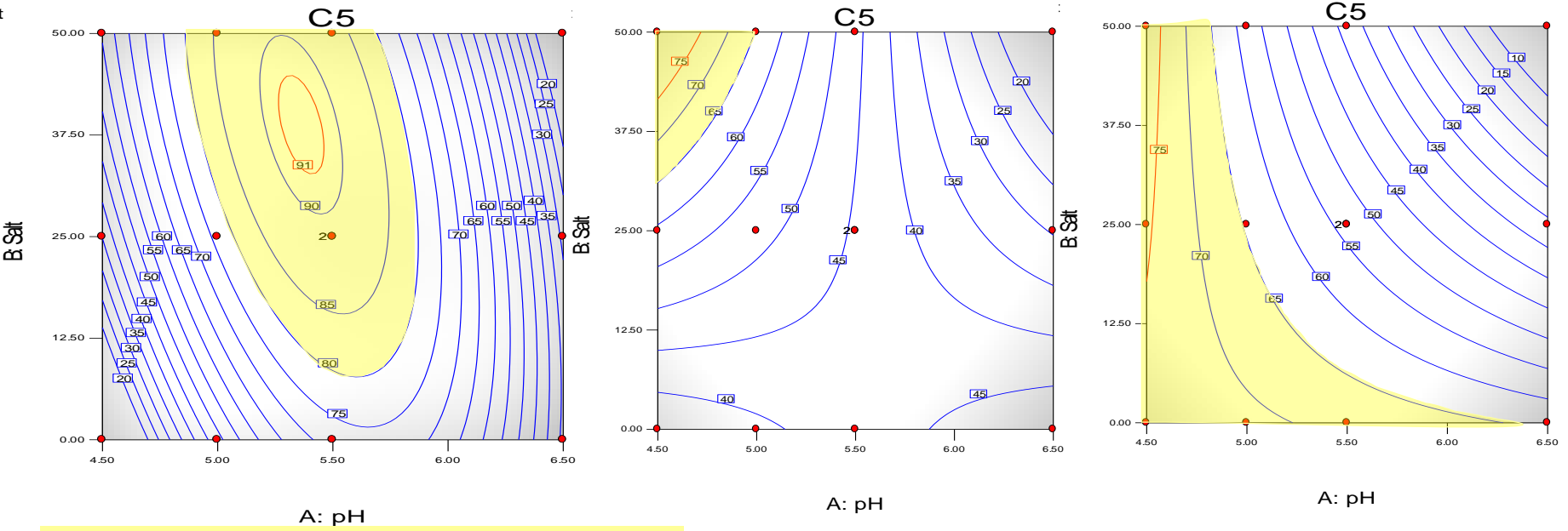

Shaded Area = 10mg/ml within max IgG C5

CEX IgG binding performance is highly variable despite all chemistries being SO_3 derived

Optimizing CEX Conditions

Binding conditions

- **pH**
 - Use a binding buffer pH 1 to 3 units below the isoelectric point of the target molecule. Dynamic binding capacity typically increases as the loading pH decreases
- **Buffer system**
 - MES, acetate, phosphate, citrate, and citrate-acetate are often used
 - Consider molecule stability, binding optimization, and the ability of the buffer to control pH in the desired operating range
- **Conductivity**
 - The load conductivity should be between 2 and 15 mS/cm
 - The optimum buffer condition depends on the target molecule and buffer pH
- **Flow rate**
 - The target operating flow rate is flexible



Optimizing CEX Conditions

Elution conditions

- **Salt gradient**

- Elution optimization should begin with a gradient elution
- Run a 20 CV gradient from low salt to approximately 500 mM to 1M sodium chloride
- Once elution performance is understood, a step elution can be implemented

- **pH**

- Start with the elution pH matched to the binding pH
- Then optimize the pH of the elution buffer as these can differ

- **Dynamic binding capacity**





- Separation as a function of DBC should be assessed
- The max DBC at which a given separation can be obtained depends on a number of factors (i.e. sample solubility, column selectivity, buffer pH, and buffer conductivity)

- **Bed height**

- Initial screening can be run with shorter bed heights, but development should be conducted at the final desired bed height


Overview

- CEX introduction: POROS® HS and XS
- Applications and Characterization Data
 - Applications
 - Viral clearance
 - Pressure vs flow curves
 - Particle size analysis
- How to optimize
- Case Studies:
 - POROS® HS process redesign
 - POROS® HS in HTP screening and scale up study
 - Rethinking polish applications: CEX in flow through/overload mode
- Summary

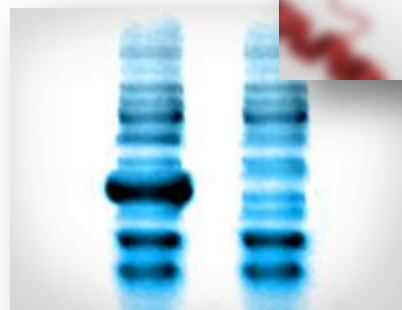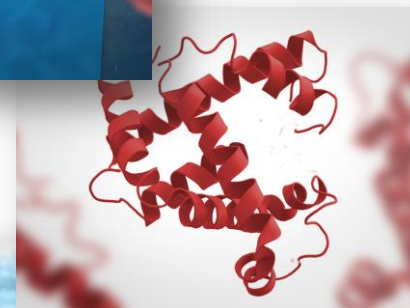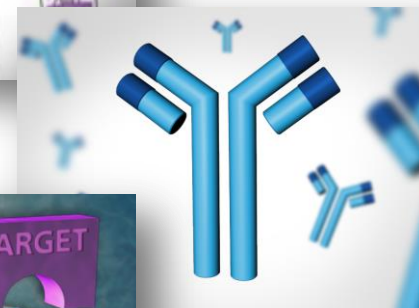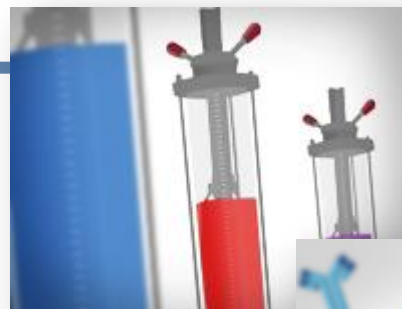

POROS[®] HS Case Study: Process Redesign

- Goal: To improve process yield and product purity within the manufacturing confines of multiple facilities for Avastin, an approved, commercially-available monoclonal antibody
- Data presented by Debbie O'Connor at May 06 IBC Antibody Development and Manufacturing Conference
 - Summarized the outcome of 2 years of process re-development by Genentech
 - 3 CEX Resins Evaluated
 - SP Sepharose Fast Flow (SCX, GE)
 - CM Sepharose Fast Flow (WCX, GE)
 - POROS HS50 (SCX, AB)
 - Protein MW 149 kDa, pI 8.1
 - Each resin loaded to 40 mg protein per mL resin
- CM Seph FF used in the existing process
- Each resin was optimized individually and the three optimized chromatography unit operations were summarized and compared

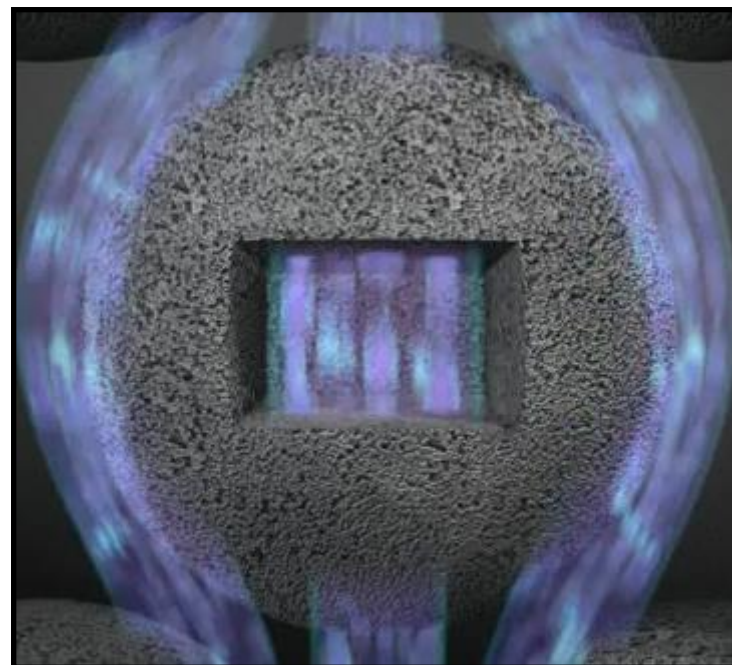


POROS[®] HS Case Study: Process Redesign

- Evaluated the following for each resin:
 - Ease/Robustness of Elution
 - Yield
 - Elution Pool Volume
 - Impurity Removal (Aggregate, HCP, DNA, Leached rPA, Retrovirus)
- POROS HS50 resulted in the simplest unit operation
 - NaCl step elution
 - Tight elution peak resulted in smallest elution peak volume which contributed to ability to scale into most stringent facility (tank size limited)
 - No temperature effects; Temp 15-25°C, no change in elution peak volume or impurity levels
 - NaCl could be used to elute the SPFF column, but [2x] required and the elution pool volume was very high

Resin	Elution Buffer	Peak Pooling
CMFF	Sodium Sulfate	Equation
SPFF	Sodium Acetate	High OD Calculation
HS	Sodium Chloride	OD

POROS[®] HS Case Study: Process Redesign

- POROS HS was chosen over SPSEphFF and CMSEphFF due to the following:
 - 11% increase in product yield
 - Improved resolution of impurities
 - HCP Removal
 - CMSEphFF << SPSEphFF < POROS HS
 - Aggregate Removal
 - SPSEphFF: 0.2-29% increase in [aggregate] across elution pool
 - POROS HS: No increase in [aggregate] across elution pool
 - Leached rPA and DNA both below LOD for POROS HS
 - Efficient binding and eluting characteristics
 - Comparatively, POROS HS could be loaded over a wider range of pH and conductivity without impacting yield or impurity removal (aggregate)
 - POROS HS presented with the same re-use, cleanability and sanitization profile as the current resin
 - Better viral clearance capability
 - 4 Log POROS HS vs. 2 log SPSEphFF
 - Resulted in the simplest and most robust step

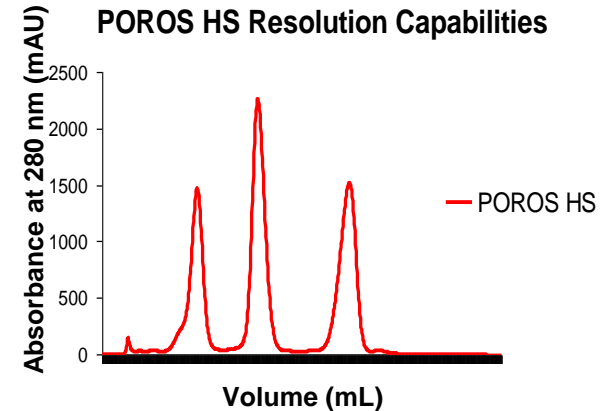

Case Study: Boehringer Ingelheim Poster Presentation at Recovery XIII Conference

- This poster was presented in June 2008 by Rathjen, Wenzel, Nothelfer, Eckermann, and Ambrosius from BI.
- RAPPTor technology (Rapid Automated Protein Purification Technology) allows for multiple binding and elution conditions to be screened and evaluated in parallel.
- Although the poster was on the applications of the screening system, the data shows the excellent performance of POROS HS compared to other resin.

BI Case Study: Purification of 2 Monoclonal Antibodies using POROS HS 50

- BI conducted an evaluation of CEX resins for recovery and Host Cell Protein removal for 2 Protein A purified mAbs
 - mAb A: human IgG1, 145 kDa, PI=7.3
 - mAb B: human IgG1, 145 kDa, PI=7.9
- 6 CEX resins were tested and 3 presented: POROS HS, Fractogel S03, and Toyopearl SP 650M
- Evaluated different pH conditions using an acetate/tris buffer system: 4.5, 5.0, 5.5, and 6.0 (Binding/eluting at same pH condition)
- Elution carried out by increasing salt concentration from 50, 100, 200 and 300 mM

BI Case Study: Purification of 2 Monoclonal Antibodies using POROS HS

- POROS HS received the highest rating for combined HCP clearance and recovery for both antibodies in the RAPPTor screening study
- Overall HS had the highest HCP removal factor of all the resins (pH 4.5: 18.9 fold HCP removal)


CEX resin	cand.	binding pH	Product recovery [%]	HCP removal factor	overall rating
Poros 50 HS	A	4.5	24	5.6	
		5.0	87	12.6	
		5.5	100	8.7	
		6.0	92	7.0	
Fractogel SO3		4.5	21	1.1	
		5.0	87	3.6	
		5.5	101	4.7	
		6.0	100	3.9	
Toyopearl SP	4.5	78	1.7		
	5.0	94	3.8		
	5.5	98	3.9		
	6.0	59	3.7		
Poros 50 HS	B	4.5	12	18.9	
		5.0	82	17.0	
		5.5	104	10.0	
		6.0	103	7.7	
Fractogel SO3		4.5	10	0.9	
		5.0	84	2.9	
		5.5	105	8.7	
		6.0	107	7.8	
Toyopearl SP	4.5	65	1.7		
	5.0	95	7.4		
	5.5	104	7.5		
	6.0	88	4.3		

Summary of CEX Screening: Product recoveries >90% and HCP removal >8 fold are highlighted in green.



BI Case Study: Purification of 2 Monoclonal Antibodies using POROS HS 50

- The RAPPTor POROS HS process was scaled up 330X by CV and yielded 93% recovery and a 7.1 fold removal of HCP showing the predictability of the model
- Column Dimensions:
1 cmD x 21 cmH, 16.5 ml
- Load, 4.6 mg mAb B/ ml resin
- Wash: 100 mM NaCl, pH 5.5
- Elution: 300 mM NaCl, pH 5.5

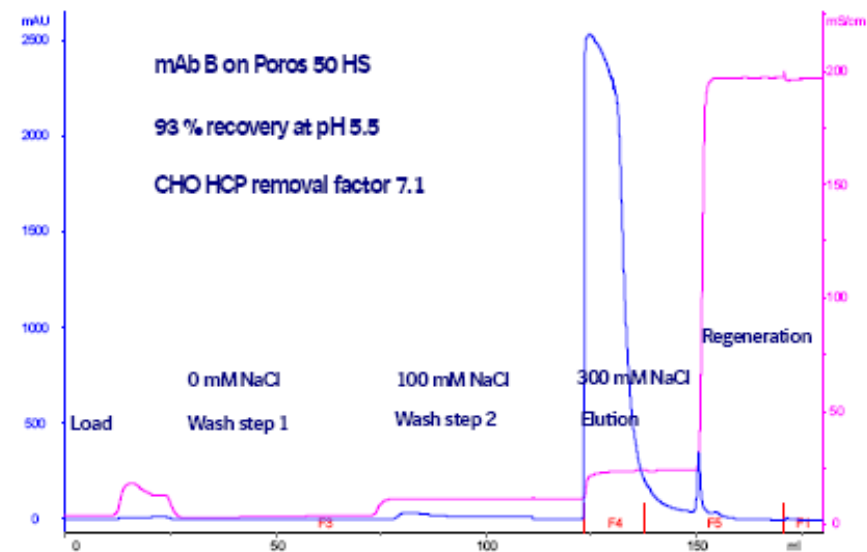

Membrane Adsorbers vs. Traditional Chromatography in Flow Through Mode

- Membrane Advantages
 - Faster flow rates
 - 95% less buffer usage
 - Disposability
 - Ease of use
- Membrane Disadvantages
 - Complicated scale-down models
 - Difficult to compare to other membranes and resin
- Chromatography Advantages
 - Reusable at commercial scale
 - Available for all scales
- Chromatography Disadvantages
 - Very large diameter columns to ↑ flow rate to ↓ bottleneck
 - Packed for throughput and not for binding capacity


So...Where does POROS® resin fit in???



POROS[®] HS Case Study: Overloaded CEX

- Goal: To evaluate CEX in overload mode to increase process efficiency and yield and reduce production costs
- Current Mab platform process: ProA > CEX (B/E) > AEX (FT)
- Although new CEX resins allow for 100 mg/ml with no loss in purity, still want to increase throughput
- Evaluated CEX in flow through type mode at binding conditions for selective displacement of impurities
 - 4 CEX products evaluated
 - POROS HS (Life)
 - SP Sepharose Fast Flow (GE)
 - Mustang S (Pall)
 - Monolith SO₃ (BIA)
 - 3 IgG1 Mabs tested with isoelectric points ranging from 8.2 to 9.2 with varying impurity amounts
 - Evaluated the following for each resin:
 - Dynamic binding capacity
 - Yield
 - Impurity Removal: HMW (dimer and higher order aggregates), CHOP, DNA, Leached rPA, Gentamicin
 - Optimized loading conditions: pH 5.5, conductivity 5 mS/cm, 600 g/L

POROS[®] Chromatography resin Impurity Breakthrough Curves

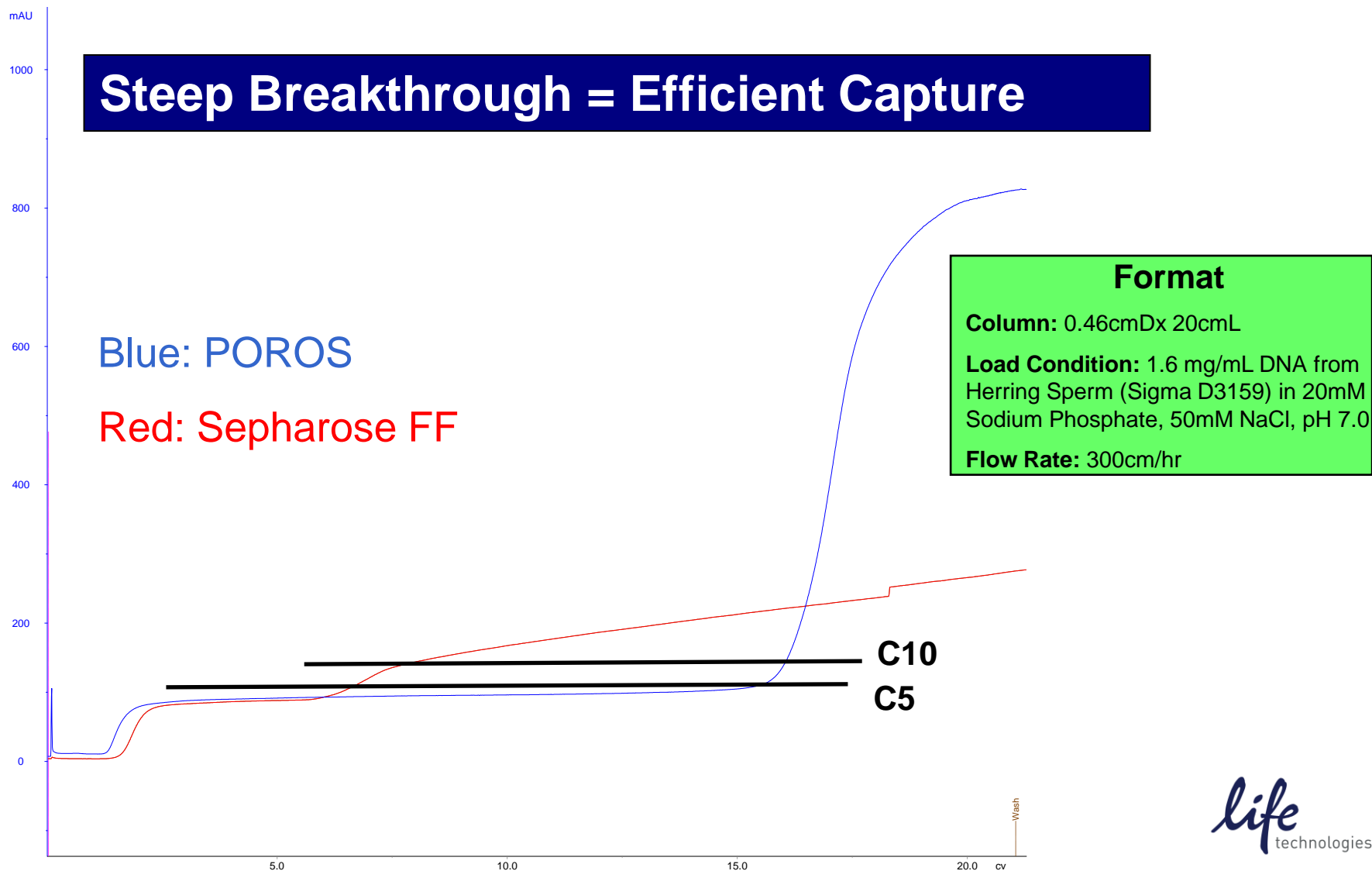


life
technologies

POROS[®] Chromatography resin: Impurity Breakthrough Curves

— POROS[®] = Dotted Lines Seph FF = Solid Lines

Steep Breakthrough = Efficient Capture

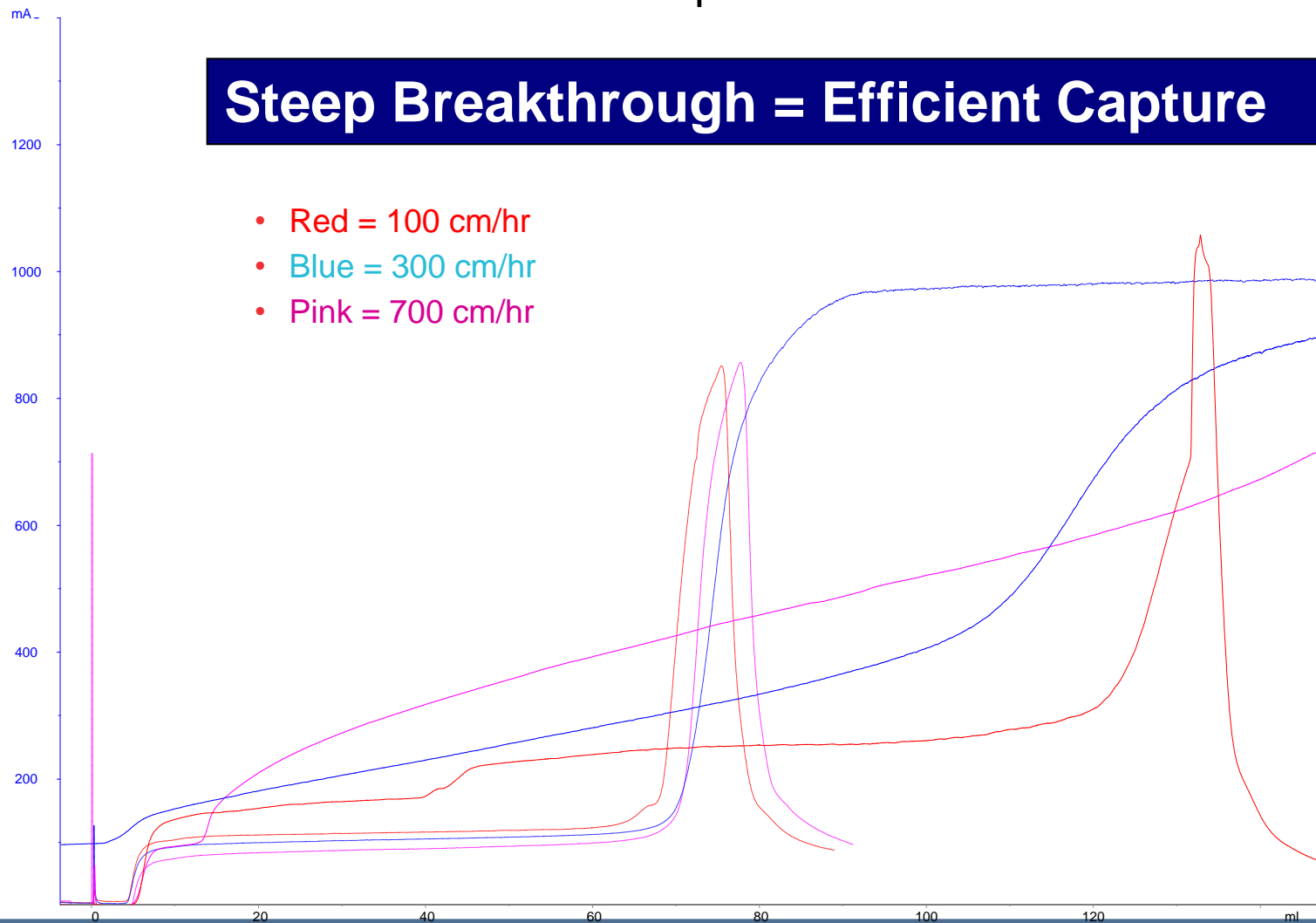

life
technologies

POROS[®] HS Case Study: Overloaded CEX

- POROS[®] HS50 resulted in the best performance
 - >90% yield when loaded at a capacity of 600 g/L
 - The best impurity removal of all CEX products especially HMW (dimers and higher order aggregates) and DNA
 - Maintained performance up to 300 cm/hr whereas SPFF decreased in performance even as low as 200 cm/hr (higher flow rates than 300 cm/hr not tested)
 - Leached ProA and Gentamicin reduced to below limit of detection
- Data based on two sources:
 - Publication: Liu, Hui F, et al, “Exploration of overloaded cation exchange chromatography for monoclonal antibody purification”, Journal of Chromatography A, V1218, 2011, Pages 6943-6952
 - Data presented by Hui Liu Ph.D., Genentech Process Development Scientist at IBC Antibody Development & Production Conference, Bellevue, Washington, March 17, 2011



Comparison of Overload vs B/E CEX Process

(13.5 kg, 2700L Process, Life Tech Model Process)

	POROS HS in Overload Format	Conventional Soft Gel Resin in B/E	High Capacity CEX in B/E	% Difference w/ POROS HS	% Difference w/ High Capacity
Column/Membrane	60 x 8 cm	140 x 22 cm	100 x 17 cm		
Column/Membrane Volume	23	339	134	93%	61%
Load Capacity (mg of protein/ml of resin)	600	40	100		
Linear Flow Rate (cm/hr or	300	100	300		
Volumetric Flow Rate (L/min)	14	26	39		
Product Load Process Time	3.2	1.7	1.2		
Total Process Time (hr)*	4.3	6.6	2.8	35%	57%
Buffer Volume (L)	542	8129	3204	93%	61%

- Both POROS HS in overload mode and High Capacity resins, such as POROS XS, in B/E mode could result in significantly smaller columns than the conventional B/E process resulting in smaller columns and equipment and less buffer usage
- Both POROS resins allow for faster flow rates resulting in shorter processing times
- Both POROS processes allow for more efficient and/or cost effective approaches to traditional B/E CEX processes
- POROS XS may offer more process flexibility to the traditional process due to better salt tolerance

Cost Comparison of an Overload vs B/E CEX Process (13.5 kg, 2700L Process, Life Tech Model Process)

	POROS HS in Overload Format	Conventional Soft Gel Resin in B/E	High Capacity CEX in B/E	Fold Savings w/ POROS HS	Fold Savings w/High Capacity CEX
Buffer Cost (\$)	1627	24386	9612	15	3
Process Labor Costs (\$)	4876	5572	4438	1	1
Cost of Resin/Membrane (\$)	45,200	338,700	393,825	7	1
Total Cost of Processing/Cycle (\$)					
1 Cycle	51,703	368,658	407,875	7	1
5 Cycles	15,543	97,698	92,815	6	1
10 Cycles	11,023	63,828	53,433	6	1
50 Cycles	7,407	36,732	21,927	5	2
100 Cycles	6,955	33,345	17,989	5	2

- Both POROS processes allow for a significant decrease in buffer costs compared to the traditional CEX process
- The POROS HS overload process could result in major cost savings due to the 93% decrease in column size and buffer usage

CEX in Flow-Through Applications - Conclusions





- Flow-Through Applications

- High impurity binding capacity and clearance over a wide range of process conditions
- Reduce resin volume, equipment sizes, processing time and process complexity
- Membrane performance in a bead
 - > Flow rate-independent performance
 - > Improved throughput potential
 - > Smaller column sizes possible
- Competitive with performance of membrane adsorbers
 - > Reusable at commercial scale
 - > Easily scalable
 - Design for commercial scale from beginning to reduce redevelopment costs
 - > GoPure™ Pre-Packed Columns for convenience and disposability



Overview

- CEX introduction: POROS® HS and XS
- Applications and Characterization Data
 - Applications
 - Viral clearance
 - Pressure vs flow curves
 - Particle size analysis
- How to optimize
- Case Studies:
 - POROS® HS process redesign
 - POROS® HS in HTP screening and scale up study
 - Rethinking polish applications: CEX in flow through/overload mode
- Summary

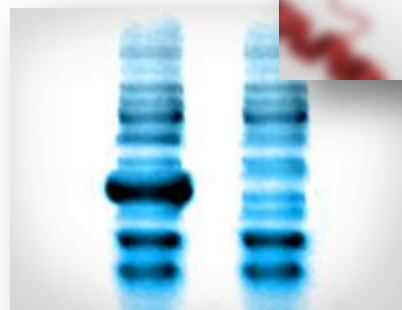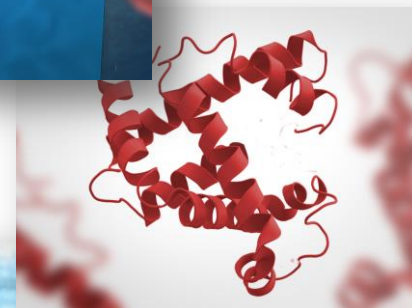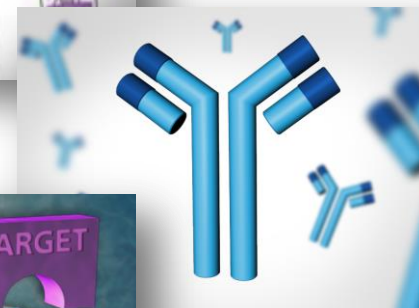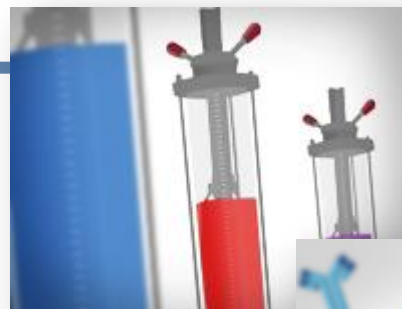

How does POROS[®] resin improve process performance and productivity?

PRODUCT ATTRIBUTE

- High Capacity



- High Resolution



- Salt Tolerance



- Efficient Chromatography



PROCESS BENEFIT

- Reduced Column Size

- Better Impurity Clearance

- Higher Product Yield

- Step versus Gradient Elution

- Reduced or Eliminated Dilution


- Elimination of Unit Operations (TFF)

- Improved Process Flow

- Reduced Pool Volume (Tank Flexibility)


POROS® XS in the Industry

- POROS® XS has been incorporated into 7 Phase 3 processes
- POROS® XS has become a platform CEX polish resin at large biopharma companies in the US and Europe
- POROS® XS is enabling non-protein A capture chromatography of monoclonal antibodies
- Customer Quotes on POROS® XS...
 - POROS® XS “provides performance consistent with your literature, POROS® XS has a higher salt tolerance, elutes at a higher salt concentration than all resins tested, and has a tight, symmetrical peak shape. Capacity looks very good”.
 - POROS® XS “has enabled a 4x reduction in process time due to capacity and throughput capability”
 - “We were able to remove 2 UF/DF steps due to POROS® XS”
 - “We incorporated POROS® XS into a troublesome biosimilar process in the course of 1 week due to the capacity advantage offered”

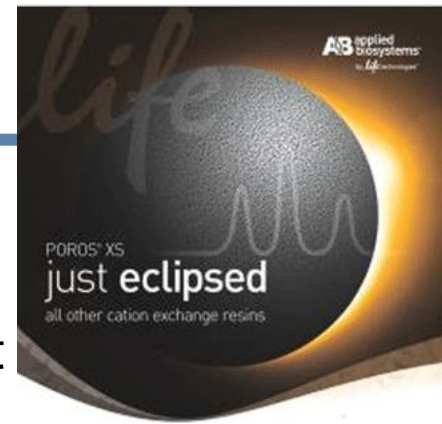

Summary


- POROS® XS sets the bar in fulfilling the most important demands in today's high titer downstream purification processes
- POROS® XS provides the combination of high capacity, high resolution, and high salt tolerance that surpasses all other cation exchange resins
- POROS® XS provides key downstream processing benefits including improved yield and throughput, with added process flexibility

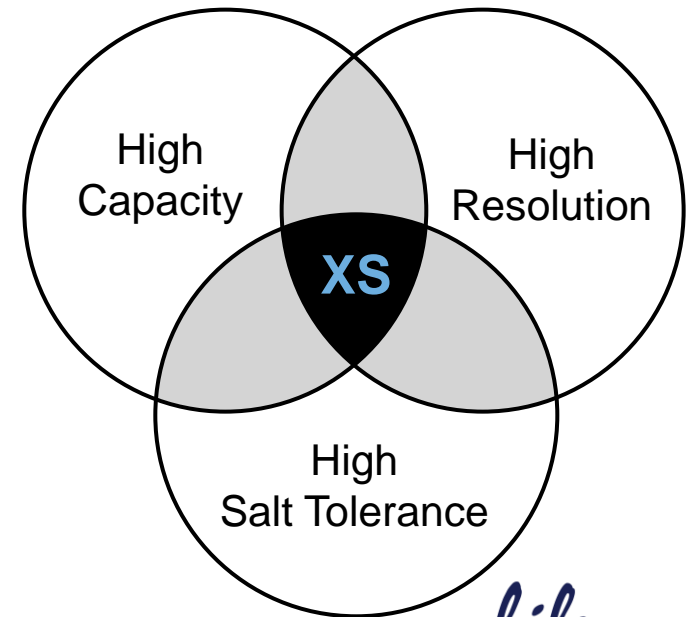

life
technologies

CEX for Intermediate Polish Chromatography

POROS® XS resin is setting standards in resolution and capacity

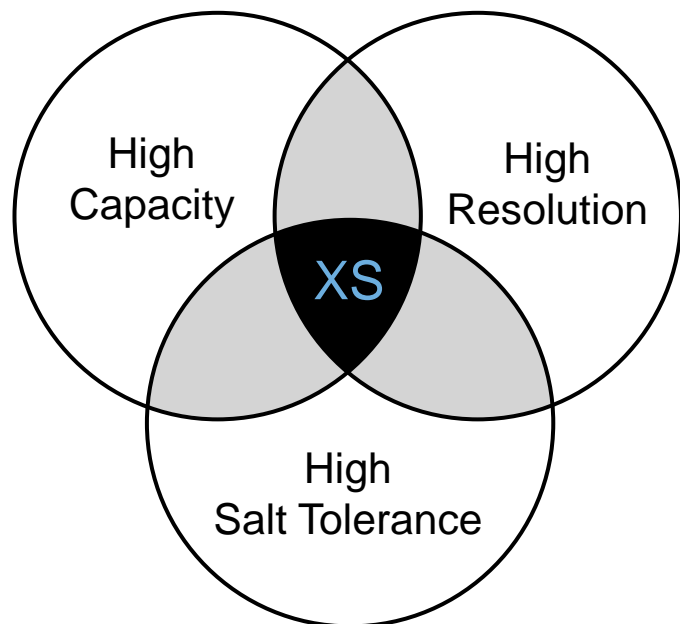

- **High capacity and resolution needed**
efficient removal of aggregate, host cell proteins, and other impurities
- **Superior resolution**
independent of load concentration or flow rate
- **More process flexibility**
over a broad range of process conditions
- **Provides solutions to today's industry demands**
POROS® XS successfully combines superior capacity, resolution and salt tolerance improving yield, throughput and increased flexibility

Features	Benefits	Products and Services
<ul style="list-style-type: none"> • High 5% breakthrough capacity >100 mg/mL over a range of process conditions • Superior resolution even with varying linear flow rates and load concentrations • Incompressible 50 µm bead with robust physical stability • Validated manufacturing process • Worldwide technical support 	<ul style="list-style-type: none"> • Improved yield • More process flexibility • Reduced column sizes • Decreased tank sizes • Improved throughput • Reduced number of unit operations • Enables CEX capture chromatography, eliminating need for rPA affinity 	<ul style="list-style-type: none"> • Bulk Resin • GoPure™ Pre-Packed Chromatography Columns