



# **5-IN-1 ELECTROPORATION RF LED FACIAL BEAUTY DEVICE**

ESRFLDFCDVA



# SAFETY & WARNINGS

Please read these safety precautions carefully and only use this device according to the instructions outlined in this guide. To make sure you can use the device correctly and avoid harm or injury to others, please follow these precautions.

## WARNING:

- Failure to follow the directions for use as they are written in this guide may lead to the device being damaged.
- If users are suffering from physical, sensory and psychiatric illness or undergoing dermatological treatments, then please use the product when you are monitored by a skin specialist, dermatologist, or your GP.
- If the charging cable is damaged, replace it immediately. Do not attempt to repair.
- Do not allow the charging cable to become damaged. This includes, but is not limited to, allowing it to bend, twist, stretch, have heavy things placed on it, or be placed in a high-temperature location.
- Do not keep the device around high temperature items such as heaters or open flames.
- Do not submerge the item in water or attempt to wash it by submerging. This device is not completely waterproof.
- This device is not to be used by children.

## ATTENTION:



- This product is only for daily use. Do not disassemble or attempt to repair it.
- For hygiene, do not allow several people to share this device. It is designed for use with one person.
- If your skin becomes irritated or hurt, cease use immediately and see your doctor.
- Do not use or store the device and its accessories in damp places (bathrooms, laundries, near showers, etc).
- If the device is damaged or broken, cease use immediately. Do not use if malfunctioning.



Do not touch the device if it is in water.

Disconnect the power to prevent from electric shock.

# OVERVIEW

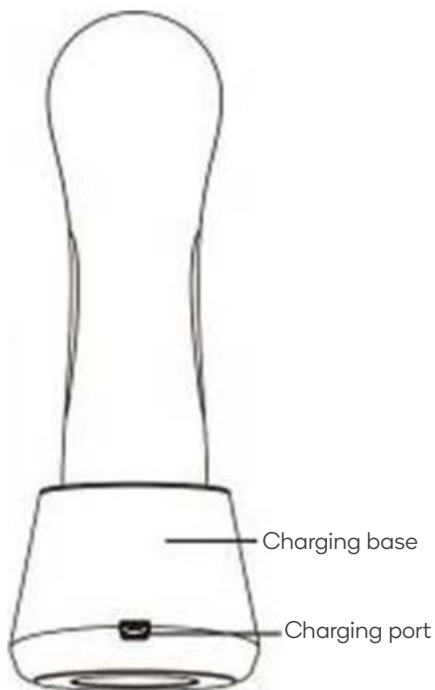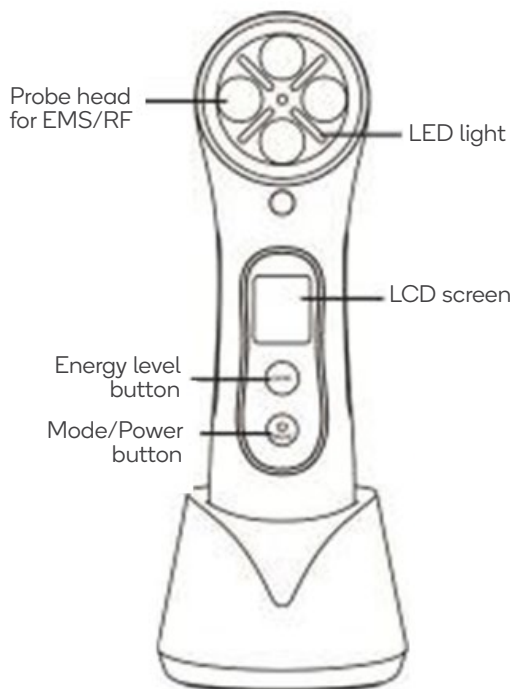

# COMPONENTS

- 1 x Facial Beauty Device
- 1 x Charging Base
- 1 x USB Charging Cable
- 1 x User Guide

# OPERATION

## Before Using

Ensure the device is charged for at least 2 hours prior to first use.

## Operating the Device






Press and hold the Mode/Power button  for 3 seconds to turn on the device.

The LCD screen will turn on and automatically enter Clean mode at level 1.

If no other mode is selected, the device will automatically switch off in 5 minutes, as this is the default Clean mode length. Other modes have longer default times.

## Energy Level Adjustment

In any mode, you can press the Level button  to cycle through the options. Levels are ordered from 1-5 and displayed by the  bars on the LCD screen.

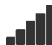LCD Display	Energy Level
	Energy Level: 1
	Energy Level: 2
	Energy Level: 3
	Energy Level: 4
	Energy Level: 5

## Operation Method

Operate the device across your face, starting from the chin and jawline, and slowly moving up to your forehead. In each area, start in the middle and move outwards in an upward motion. Ensure the skin is moisturised for better gliding movements.

Always use the product in conjunction with facial essence or moisturiser (sold separately). This device achieves best results with water-based gel moisturisers. Apply essence after use.

Recommended use is two to three times a week.

# MODES

Press Mode/Power button  to select the mode.

The selected mode will be identified by the ✓ icon on the LCD display.

The five modes use RF (Radio Frequency), EMS (Micro current), Acoustic Vibration, and LED lighting technology together:

Mode	Description/Function	Default Time
<b>Clean Mode</b> Green LED + Sonic vibration	Clean up harmful substances in the skin, reduce pigmentation, improve skin tone, accelerate blood circulation and skin metabolism.	5 minutes
<b>Lift Mode</b> Yellow LED + EMS + RF	Restore skin flexibility, improve wrinkles, massage the muscles and stimulate the lymphatic system, redefines the jawline.	10 minutes
<b>Lead in Mode</b> Pink LED + EMS + RF	Create space between cells in a short time and effectively assist in essence and nutrient absorption.	10 minutes
<b>Anti-aging Mode</b> Red LED + Sonic Vibration + RF	Promote skin collagen regeneration, enhance skin vitality, tighten skin and improve skin tone.	10 minutes
<b>Remove Acne Mode</b> Blue LED	Promote blood circulation, anti-inflammatory, sterilisation, effectively remove the skin acne and shrink pores.	5 minutes

Once you select a mode, the LCD screen will display a timer counting down the mode's length. Once the time is complete, the device will automatically turn off.

To manually turn the device off, long press Mode/Power button  for 3 seconds.

## **Function overview:**

### **RF (Radio Frequency) Function**

RF is a high-frequency thermal penetration therapy. It transmits a kind of electromagnetic energy caused by high frequency AC change to the subcutaneous tissue. It can also stimulate the dermis collagen regeneration, increasing the thickness and the density of the dermis. It will lift and fill wrinkles, eliminate the scars and restore skin elasticity and lustre.

### **EMS (Micro current) Function**

It is a special technology of electrical perforation method and it can in a short time create space between the cells so that the liquid components of essence can be directly delivered to the skin surface and dermal layer. It ensures effective absorption of moisturiser and facial essence.

### **Acoustic Vibration Function**

It can massage and smooth muscles, promote blood circulation, speed up metabolism and tighten the skin.

## **LED Functions**

### **Green light - Wavelength: 520nm (35nm)**

Reduce pigmentation, discolouration, freckles, sunspots and age spots. Improve and brighten skin tone, continuous use can prevent development of pigments.

### **Yellow light - Wavelength: 590nm (35nm)**

Stimulate the lymph and nerve system which regulate muscle elasticity. Regular use of yellow light can assist in skin tightening and better toning.

### **Pink light - Wavelength: 700nm (35nm)**

Increase blood circulation and improves metabolism, assist in effective absorption of skin care products.

### **Red light - Wavelength: 620nm (35nm)**

Promote collagen regeneration, smooth wrinkles and enhance skin vitality.

### **Blue light - Wavelength: 465nm (35nm)**

Reduce bacteria that causes skin inflammation, reduce acne and redness.

# CHARGING

## Battery status

The battery status is shown on the LCD screen with three grids in the battery icon.

LCD status	Battery remaining
1 Grid	0-30%
2 Grid	30%-70%
3 Grid	70%-100%

## Charging

- When the battery is low, the battery icon will begin flashing. We recommend charging to 100% (full three grids, no longer flashing) before further use.
- Connect the USB charging cable to the Charging Base. Connect the cable to any appropriate USB adapter (5V, 500mA) or USB port.
- The battery icon flashes during charging.
- When the charge is complete, the battery icon will stop flashing and remain solid.



# CLEANING & CARE

- Do not submerge in water. Do not use water to clean this device directly.
- Do not clean while connected directly to USB charger.
- After operation, use paper, towels or dry soft cloth to clean the head of the device.
- When the dirt is difficult to remove, try using a slightly damp cloth to wipe off.
- If not using the device for a long time, it is recommended that you charge the device regularly and store it in a dry, cool place.
- The device is not completely waterproof; please do not use water to clean the instrument to avoid damage.

# TROUBLESHOOTING

## Common Problems and Solutions

Issue	Troubleshooting
Not feeling the device on face	Ensure the device is held with the head parallel to your face.
	Adjust the energy power.
The device is not turning on	Insufficient battery charge. Charge the device.
The energy level is not available	Internal parts of the device damaged.
	Button stuck.

If you are still experiencing issues, or if your device is not operating as it should be, stop using and contact [help.kogan.com](https://help.kogan.com) for assistance.

# SPECIFICATIONS

<b>Charging Voltage</b>	3.7V, 100mA
<b>Charging Hours</b>	2 hours
<b>Battery Capacity</b>	3.7V, 830mAh
<b>Vibration Frequency</b>	8600 ( $\pm 10$ ) rpm
<b>Net Weight</b>	135g
<b>Charging Base</b>	DC 5V 500mA



### **Need more information?**

We hope that this user guide has given you the assistance needed for a simple set-up. For the most up-to-date guide for your product, as well as any additional assistance you may require, head online to **[help.kogan.com](https://help.kogan.com)**

kogan