

November 2014

Investigator[®] DIPplex Handbook

For multiplex amplification of 30
deletion/insertion polymorphisms, plus
Amelogenin



Sample & Assay Technologies

QIAGEN Sample and Assay Technologies

QIAGEN is the leading provider of innovative sample and assay technologies, enabling the isolation and detection of contents of any biological sample. Our advanced, high-quality products and services ensure success from sample to result.

QIAGEN sets standards in:

- Purification of DNA, RNA, and proteins
- Nucleic acid and protein assays
- microRNA research and RNAi
- Automation of sample and assay technologies

Our mission is to enable you to achieve outstanding success and breakthroughs. For more information, visit www.qiagen.com.

Contents

Kit Contents	4
Storage	4
Intended Use	4
Safety Information	5
Introduction	6
Equipment and Reagents to Be Supplied by User	10
Protocols	
■ PCR Amplification	11
■ Electrophoresis Using the ABI PRISM 310 Genetic Analyzer	14
■ Electrophoresis Using the ABI PRISM 3100-<i>Avant</i>/3100 Genetic Analyzer	21
■ Electrophoresis Using the Applied Biosystems 3130/3130xl Genetic Analyzer	29
■ Electrophoresis Using the Applied Biosystems 3500/3500xL Genetic Analyzer	40
■ Analysis	52
Interpretation of Results	59
Troubleshooting Guide	60
References	63
Ordering Information	64

Kit Contents

Investigator DIPplex Kit	(100)
Catalog no.	384015
Number of 25 µl reactions	100
Reaction Mix A	500 µl
Primer Mix DIPplex	500 µl
Multi Taq2 DNA polymerase	150 U
Nuclease-free water	2 x 1.9 ml
Control DNA 9948	200 µl
DNA size standard 550 (BTO)	50 µl
Allelic ladder DIPplex	25 µl
Quick-Start Protocol	1

Storage

All components of the Investigator DIPplex Kit should be stored at –30 to –15°C. Avoid repeated thawing and freezing. Primer mix and allelic ladder must be stored protected from the light. DNA samples and post-PCR reagents (allelic ladder and DNA size standard) should be stored separately from the PCR reagents. Under these conditions, the components are stable until the expiration date indicated on the kit.

Intended Use

The Investigator DIPplex Kit is intended for molecular biology applications in human identity testing. This product is not intended for the diagnosis, prevention, or treatment of a disease.

USA and Canada: not for use in, or in preparation for, legal proceedings, including parentage determination.

All due care and attention should be exercised in the handling of the products. We recommend all users of QIAGEN® products to adhere to the NIH guidelines that have been developed for recombinant DNA experiments, or to other applicable guidelines.

Safety Information

When working with chemicals, always wear a suitable lab coat, disposable gloves, and protective goggles. For more information, please consult the appropriate safety data sheets (SDSs). These are available online in convenient and compact PDF format at www.qiagen.com/safety where you can find, view, and print the SDS for each QIAGEN kit and kit component.

Introduction

Human identification is commonly based on the analysis of short tandem repeats (STRs) or single nucleotide polymorphisms (SNPs), depending on the demands of an examination or on the sample quality. Analyzing biallelic deletion/insertion polymorphisms (DIPs) combines the advantages of both STR and SNP analysis. DIPs are particularly useful for identification in anthropological or population genetic applications. In addition, DIP analysis confers further benefits, such as the absence of stutter peaks, enabling mixed stain analysis.

The Investigator DIPplex Kit was developed for multiplex amplification of 30 biallelic DIPs plus Amelogenin. The 30 DIPs are distributed over 19 autosomes that are at least 10 Mbp away from any commercially available STR and SNP marker (Table 2, page 8). Therefore, in combination with standard markers, these DIP markers improve discriminatory power.

The primers are fluorescence-labeled with one of the following dyes:

- 6-FAM™: Amelogenin, HLD77, HLD45, HLD131, HLD70, HLD6, HLD111, HLD58, HLD56
- BTG: HLD118, HLD92, HLD93, HLD99, HLD88, HLD101, HLD67
- BTY: HLD83, HLD114, HLD48, HLD124, HLD122, HLD125, HLD64, HLD81
- BTR: HLD136, HLD133, HLD97, HLD40, HLD128, HLD39, HLD84

We recommend using DIPSorter Freeware for easy sorting of genotypes and interpretation of results.

Table 1. Discrimination power of DIPs, STRs, and SNPs

	Loci	CPE/Trio*	CPI†	Population
Investigator DIPplex	30 DIPs	0.9980	2.83×10^{-13}	German Caucasian
AmpF ℓ STR® Minifiler™	8 STRs	0.99976	8.21×10^{-11}	US Caucasian
AmpF ℓ STR SEfiler Plus™	11 STRs	0.999998	7.46×10^{-14}	US Caucasian
Powerplex® 16	15 STRs	0.9999994	5.46×10^{-18}	US Caucasian
Sanchez et al. 2006 (1)	52 SNPs	0.9998	5.00×10^{-21}	European Caucasian

* Combined probability of paternity exclusion.

† Combined probability of identity.

The optimal amount of DNA under standard conditions is 0.2–0.5 ng. Internal validations demonstrated reliable results with <0.1 ng DNA.

The Investigator DIPplex Kit was validated using the GeneAmp® PCR System 9700 (with Gold-plated 96-Well Silver Block), ABI PRISM® 310, ABI PRISM 3100, Applied Biosystems 3130, and Applied Biosystems 3500 Genetic Analyzers.

Table 2 (page 8) shows the DIP loci of the Investigator DIPplex, their chromosomal localization, GenBank®/SNP ID number, motif, and respective reference allele.

Table 2. Locus-specific information of the Investigator DIPplex Kit

DIP locus	Chromosome localization	GenBank accession/ SNP ID	Motif (+DIP)	Reference allele
DIPplex Blue				
Amelogenin X	Xp22.1-22.3	M55418		X
Amelogenin Y	Yp11.2	M55419		Y
HLD77	7q31.1	rs1611048	TAAG	+DIP
HLD45	2q31.1	rs2307959	CACG	–DIP
HLD131	7q36.2	rs1611001	TGGGCTTATT	+DIP
HLD70	6q16.1	rs2307652	AGCA	–DIP
HLD6	16q13	rs1610905	GCAGGACTGG GCACC	–DIP
HLD111	17p11.2	rs1305047	CACA	–DIP
HLD58	5q14.1	rs1610937	AGGA	+DIP
HLD56	4q25	rs2308292	TAAGT	+DIP
DIPplex Green				
HLD118	20p11.1	rs16438	CCCCA	–DIP
HLD92	11q22.2	rs201771066	GTTT	–DIP
HLD93	12q22	rs150042219	ACTTT	–DIP
HLD99	14q23.1	rs2308163	TGAT	–DIP
HLD88	9q22.32	rs8190570	CCACAAAGA	+DIP
HLD101	15q26.1	rs2307433	GTAG	–DIP
HLD67	5q33.2	rs1305056	CTACTGAC	–DIP

Table continued on next page; HLD: Human locus DIP; –DIP: Deletion; +DIP: Insertion.

Table continued from previous page.

DIP locus	Chromosome localization	GenBank accession/ SNP ID	Motif (+DIP)	Reference allele
DIPplex Yellow				
HLD83	8p22	rs2308072	AAGG	–DIP
HLD114	17p13.3	rs2307581	TCCTATTCTACT CTGAAT	–DIP
HLD48	2q11.2	rs28369942	GACTT	–DIP
HLD124	22q12.3	rs6481	GTGGA	–DIP
HLD122	21q22.11	rs8178524	GAAGTCTGAGG	–DIP
HLD125	22q11.23	rs16388	ATTGCC	–DIP
HLD64	5q12.3	rs397832668	GACAAA	+DIP
HLD81	7q21.3	rs17879936	GTAAGCATTGT	–DIP
DIPplex Red				
HLD136	22q13.1	rs16363	TGTTT	–DIP
HLD133	3p22.1	rs2067235	CAACCTGGATT	
HLD97	13q12.3	rs17238892	AGAGAAAGCTG AAG	–DIP
HLD40	1p32.3	rs146044344	GGGACAGGTGG CCACTAGGAGA	+DIP
HLD128	1q31.3	rs2307924	ATTAAATA	–DIP
HLD39	1p22.1	rs17878444	CCTAAACAAAAA TGGGAT	–DIP
HLD84	8q24.12	rs3081400	CTTTC	–DIP

HLD: Human locus DIP; –DIP: Deletion; +DIP: Insertion.

Equipment and Reagents to Be Supplied by User

When working with chemicals, always wear a suitable lab coat, disposable gloves, and protective goggles. For more information, consult the appropriate safety data sheets (SDSs), available from the product supplier.

- Hi-Di™ Formamide, 25 ml (Applied Biosystems®, cat. no. 4311320)
- Matrix Standards BT5 for single-capillary instruments, e.g., ABI PRISM 310 Genetic Analyzer (QIAGEN, cat. no. 386113)
- Matrix Standards BT5 for multi-capillary instruments, e.g., ABI PRISM 3100 and Applied Biosystems 3130 and 3500™ Genetic Analyzers (QIAGEN, cat. nos. 386123 or 386125)
- Pipets and pipet tips
- One of the following DNA analyzers:
 - ABI PRISM 310 Genetic Analyzer
 - ABI PRISM 3100-Avant™/3100 Genetic Analyzer
 - ABI PRISM 3130/3130xl Genetic Analyzer
 - Applied Biosystems 3500/3500xL Genetic Analyzer
- One of the following PCR thermal cyclers:
 - QIAGEN Rotor-Gene® Q
 - GeneAmp PCR System 9700
 - Bio-Rad PTC-200
 - Techne TC-512
 - Biometra T1
 - Eppendorf® Mastercycler® ep
- PCR tubes or plates

Validity analysis software for human identification products

Investigator Human Identification PCR Kits require calibration with an allelic ladder. Therefore, the software used must be compatible with human identification (HID) products for forensic applications. We recommend GeneMapper® ID, GeneMapper ID-X, or Genotyper® Software. The Investigator Template Files facilitate data analysis and are valid with the software mentioned above.

Protocol: PCR Amplification

This protocol is for PCR amplification of DIP loci from human DNA samples using the Investigator DIPplex Kit.

Important points before starting

- Set up all reaction mixtures in an area separate from that used for DNA isolation and PCR product analysis (post-PCR).
- Use disposable tips containing hydrophobic filters to minimize cross-contamination.

Things to do before starting

- Before opening the tubes with PCR components, vortex, and then centrifuge briefly to collect contents at the bottom of the tubes.

Procedure

1. Thaw PCR components and template nucleic acid.

Mix thoroughly before use.

2. Prepare a master mix according to Table 3, page 12.

The master mix contains all of the components needed for PCR except the template (sample) DNA and nuclease-free water.

Prepare a volume of master mix 10% greater than that required for the total number of PCR assays to be performed. This should include positive and negative control reactions.

3. Mix the master mix thoroughly and dispense appropriate volumes into PCR tubes or the wells of a PCR plate.

4. Add template DNA and nuclease-free water to the master mix, to give a final sample volume of 25 μ l.

5. Prepare positive and negative controls.

Positive control: Use 5 μ l of the Control DNA.

Negative control: Use nuclease-free water instead of template DNA in the reaction.

Table 3. Reaction setup

Component	Volume per reaction
Reaction Mix A*	5.0 μ l
Primer Mix	5.0 μ l
Multi Taq2 DNA Polymerase	0.6 μ l
Nuclease-free water (added in step 4)	Variable
Template DNA (added in step 4)	Variable
Total volume	25 μl

* Contains dNTP mix, MgCl₂, and bovine serum albumin (BSA).

6. Program the thermal cycler according to the manufacturer's instructions, using the conditions outlined in Table 4.

For stains containing small amounts of DNA (<100 pg/25 μ l reaction), we recommend using the cycling conditions outlined in Table 5 (page 13).

Note: If using the GeneAmp PCR System 9700 with an Aluminum Block, use "Std Mode", or with a Silver 96-Well Block or Gold-plated Silver 96-Well Block, use "Max Mode". Do not use "9600 Emulation Mode".

Table 4. Standard cycling protocol, recommended for all DNA samples

Temperature	Time	Number of cycles
94°C*	4 min	–
94°C	30 s	30 cycles
61°C	120 s	
72°C	75 s	
68°C	60 min	–
10°C	∞	–

* Hot-start to activate DNA polymerase.

Table 5. Optional cycling protocol, recommended for stains containing small amounts (<100 pg) of DNA

Temperature	Time	Number of cycles
94°C*	4 min	–
94°C	30 s	32 cycles
61°C	120 s	
72°C	75 s	
68°C	60 min	–
10°C	∞	–

* Hot-start to activate DNA polymerase.

- 7. After the cycling protocol is completed, store samples at –30 to –15°C protected from the light, or proceed directly with electrophoresis.**

Protocol: Electrophoresis Using the ABI PRISM 310 Genetic Analyzer

For general instructions on instrument setup, matrix generation, and application of the GeneScan® or GeneMapper *ID* Software, refer to the *ABI PRISM 310 Genetic Analyzer User's Manual*. Electrophoresis using the GeneScan Software is described below.

The virtual filter set G5 is used for combined application of the 5 fluorescent labels 6-FAM, BTG, BTY, BTR, and BTO. This matrix standard is known as BT5.

The materials required for electrophoresis are given in Table 6.

Table 6. Materials required for electrophoresis

Material	Specifications
Capillary	47 cm/50 μ m (green)
Polymer	POP-4™ for ABI PRISM 310 Genetic Analyzer
Buffer	10x Genetic Analyzer Buffer with EDTA

Matrix generation

Before conducting DNA fragment size analysis with the filter set G5, a matrix with the 5 fluorescent labels 6-FAM, BTG, BTY, BTR, and BTO must be generated (Table 7).

Table 7. The fluorescent labels of BT5

Color	Matrix standard
Blue (B)	6-FAM
Green (G)	BTG
Yellow (Y)	BTY
Red (R)	BTR
Orange (O)	BTO

1. **Five electrophoresis runs should be conducted, one for each fluorescent label, under the same conditions as for the samples and allelic ladders of the Investigator DIPplex Kit, in order to generate suitable matrix files (Table 8).**

Table 8. Matrix setup for a single capillary instrument (ABI PRISM 310 Genetic Analyzer)

Matrix sample	Component	Volume
Matrix sample 1	Hi-Di Formamide	12 μ l
	Matrix standard 6-FAM	1 μ l
Matrix sample 2	Hi-Di Formamide	12 μ l
	Matrix standard BTG	1 μ l
Matrix sample 3	Hi-Di Formamide	12 μ l
	Matrix standard BTY	1 μ l
Matrix sample 4	Hi-Di Formamide	12 μ l
	Matrix standard BTR	1 μ l
Matrix sample 5	Hi-Di Formamide	12 μ l
	Matrix standard BTO	1 μ l

2. **Denature for 3 min at 95°C.**
3. **Snap freeze by placing the plate on ice for 3 min.**
Alternatively, the thermal cycler set to 4°C may be used to cool the plate.
4. **Load the samples on the tray.**
5. **Create a Sample Sheet and enter the sample designation. Table 9 (page 16) shows the injection list for matrix generation.**

Table 9. Injection list for matrix generation


Parameter	Settings
Module File	GS STR POP-4 (1 ml) G5
Matrix File	None
Size Standard	None*
Injection Time (s)	5
Injection Voltage (kV)	15
Run Voltage (kV)	15
Run Temperature (°C)	60
Run Time (min)	24

* Always prepare matrix standards without DNA Size Standard (BTO).

Analysis of the matrix samples

1. Run the GeneScan Software.
2. Select "New" from the File menu, and then select "Project".
3. Open the folder of the current run and select "Add Sample Files".
4. Select a matrix sample in the "Sample File" column.
5. Click "Sample" and then "Raw Data".
6. Check the matrix samples for a flat baseline. As shown in the figure (next page), there should be at least 5 peaks with peak heights of 1000–4000 RFU for each matrix sample.

Note: The optimal range is 2000–4000 RFU.

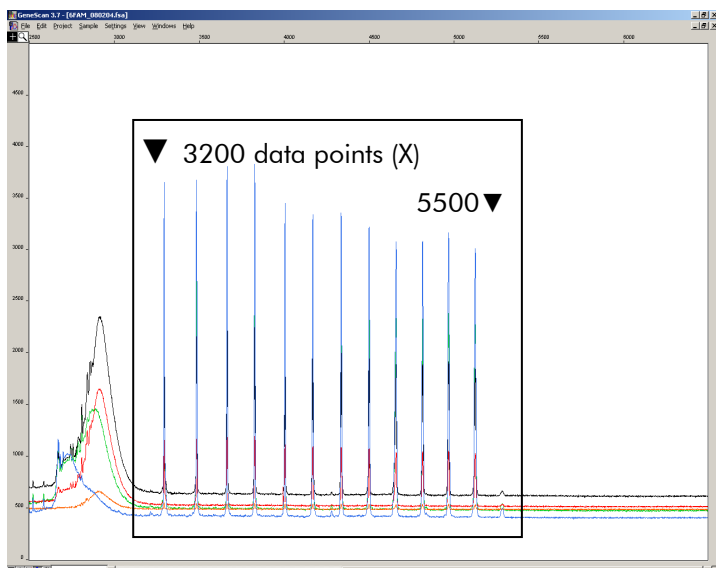


Electropherogram with raw data of the matrix standard 6-FAM.

- 7. Select an analysis range with a flat baseline and re-inject the matrix sample, if necessary.**
- 8. Record start and end values (data points) of the analysis range; e.g., start value 3200, end value 5500.**
- 9. Calculate the difference between the end and start values; e.g., $5500 - 3200 = 2300$ data points.**


Generation of a matrix

- 1. Select "New" in the "File" menu, and then select "Matrix".**
- 2. Import the matrix samples for all dyes (B, G, Y, R, and O).**
- 3. Enter a "Start At" value, e.g., 3200.**
- 4. Under "Points", enter the calculated difference between end and start values, e.g., 2300.**
- 5. Click "OK" to calculate the new matrix.**

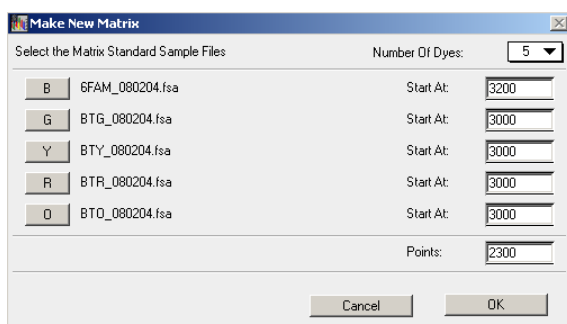

Matrix sample selection.

6. Select "Save as" in the File menu to save the new matrix in the matrix folder.

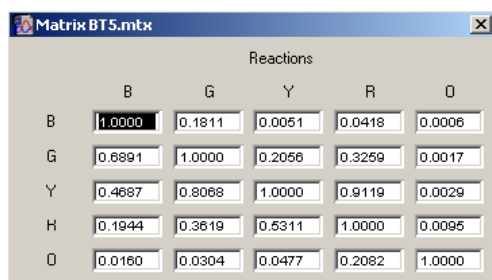

The screenshot shows a window titled "Matrix BT5.mtx" with a table of reaction values. The table has columns for dyes B, G, Y, R, and O, and rows for the same dyes. The values are as follows:

	B	G	Y	R	O
B	1.0000	0.1811	0.0051	0.0418	0.0006
G	0.6891	1.0000	0.2056	0.3259	0.0017
Y	0.4687	0.8068	1.0000	0.9119	0.0029
H	0.1944	0.3619	0.5311	1.0000	0.0095
O	0.0160	0.0304	0.0477	0.2082	1.0000

New matrix BT5.

Checking the matrix

1. To check the new matrix with current samples, select "New" in the "File" menu, and then select "Project".
2. Open the folder of the respective run and select "Add Sample Files".
3. Select the sample(s) in the "Sample File" column.
4. Click "Sample" and then "Install New Matrix" to open the matrix folder and select the new matrix.
5. Re-analyze the samples.

Note: There should be no pull-up peaks between the dye panels (B, G, Y, R, O) with the new matrix.

Sample preparation

1. Set up a mixture of formamide and DNA size standard according to Table 10.

Table 10. Setup of formamide and DNA size standard mixture

Component	Volume per sample
Hi-Di Formamide	12.0 μ l
DNA Size Standard 550 (BTO)	0.5 μ l

2. For each sample to be analyzed, aliquot 12 μ l of the mixture to a tube.
3. Add 1 μ l PCR product or allelic ladder (diluted, if necessary).
4. Denature for 3 min at 95°C.

5. Snap freeze by placing the plate on ice for 3 min.

Alternatively, the thermal cycler set to 4°C may be used to cool the plate.

6. Load the samples on the tray.

Setting up the GeneScan Software

Create a Sample Sheet and enter sample designation (Table 11).

Table 11. Injection list for ABI PRISM 310 Genetic Analyzer

Component	Settings
Module File	GS STR POP-4 (1 ml) G5
Matrix File	e.g., Matrix BT5
Size Standard	e.g., SST-BTO_60-260bp
Injection Time (s)	5*
Injection Voltage (kV)	15
Run Voltage (kV)	15
Run Temperature (°C)	60
Run Time (min)	20 [†]

* Deviating from standard settings, the injection time may range between 1 and 10 s depending on the type of sample. If samples with very high signal intensities are recorded, a shorter injection time may be selected. For samples with low DNA content, an injection time up to 10 s may be necessary.

[†] The run time for Investigator DIPplex was modified in order to be able to analyze fragments with lengths of up to 200 bp.

Analysis parameters

Table 12 lists the recommended analysis parameters.

Table 12. Recommended analysis parameters for the ABI PRISM 310 Genetic Analyzer

Parameter	Settings
Analysis Range	Start: 2000 Stop: 6000
Data Processing	Baseline: Checked Multi-component: Checked Smooth options: Light
Peak Detection	Peak Amplitude Thresholds B:* Y:* G:* R:* O:* Min. Peak Half Width: 2 pts Polynomial Degree: 3 Peak Window Size: 11 pts [†]
Size Call Range	Min: 60 Max: 550
Size Calling Method	Local Southern Method
Split Peak Correction	None

* The peak amplitude threshold (cutoff value) corresponds to the minimum peak height that will be detected by the GeneScan or GeneMapper *ID* Software. Thresholds are usually 50–200 RFU and should be determined individually by the laboratory. Recommendation: The minimal peak height should be 3 times as high as the background noise of the baseline.

[†] Only the setting for Peak Window Size is different to defaults from Applied Biosystems for HID analysis.

Note: For information on using the recommended Template Files (as analysis parameters), refer to the appropriate Investigator template files user guide (Genotyper, GeneMapper *ID*, or GeneMapper *ID-X*).

Protocol: Electrophoresis Using the ABI PRISM 3100-Avant/3100 Genetic Analyzer

For detailed instructions on instrument setup, spectral calibration, application of the ABI PRISM 3100 Data Collection Software version 1.01 or 1.1 and the GeneScan Software, refer to the *ABI PRISM 3100-Avant/3100 Genetic Analyzer User's Manual*.

The system with 4 capillaries is the ABI PRISM 3100-Avant Genetic Analyzer, and the system with 16 capillaries is the ABI PRISM 3100 Genetic Analyzer.

The virtual filter set G5 is used for combined application of the 5 fluorescent labels 6-FAM, BTG, BTY, BTR, and BTO. This matrix standard is known as BT5.

The materials required for electrophoresis are given in Table 13.

Table 13. Materials required for electrophoresis

Material	Specifications
Capillary	36 cm Capillary Array for ABI PRISM 3100-Avant/3100 Genetic Analyzer
Polymer	POP-4 Polymer for ABI PRISM 3100-Avant/3100 Genetic Analyzer
Buffer	10x Genetic Analyzer Buffer with EDTA

Spectral calibration/matrix generation

Proper spectral calibration is critical for evaluation of multicolor systems with the ABI PRISM 3100-Avant/3100 Genetic Analyzer and should be done before conducting fragment length analysis. The calibration procedure creates a matrix which is used to correct the overlapping of fluorescence emission spectra of the dyes.

Spectral calibration comprises the following steps:

- Preparation of the spectral calibration standards
- Loading the standards to the 96-well reaction plate (one sample per capillary)
- Entering the plate composition
- Performing a spectral calibration run and checking the matrix

Preparing the spectral calibration standards

Example for 4 capillaries (ABI PRISM 3100-Avant Genetic Analyzer)

1. Set up a mixture of formamide and Matrix Standard BT5 according to Table 14.

Table 14. Setup of formamide and Matrix Standard BT5 mixture for 4 capillaries

Component	Volume
Hi-Di Formamide	60 μ l
Matrix Standard BT5 multi cap.	5 μ l

2. Load 12 μ l of the mixture to a 96-well plate; e.g., positions A1–D1.
3. Denature for 3 min at 95°C.
4. Snap freeze by placing the plate on ice for 3 min.

Alternatively, the thermal cycler set to 4°C may be used to cool the plate.

Example for 16 capillaries (ABI PRISM 3100 Genetic Analyzer)

1. Set up a mixture of formamide and Matrix Standard BT5 according to Table 15.

Table 15. Setup of formamide and Matrix Standard BT5 mixture for 16 capillaries

Component	Volume
Hi-Di Formamide	204 μ l
Matrix Standard BT5 multi cap.	17 μ l

2. Load 12 μ l of the mixture to a 96-well plate; e.g., positions A1–H1 and A2–H2.
3. Denature for 3 min at 95°C.
4. Snap freeze by placing the plate on ice for 3 min.

Alternatively, the thermal cycler set to 4°C may be used to cool the plate.

Performing a spectral calibration run

The parameter file for Dye Set G5 must be modified once to achieve successful calibration with the Data Collection Software version 1.0.1 or 1.1.

Spectral parameter

1. To change settings in the parameter file, go to the following path:
D:\AppliedBio\Support Files\Data Collection
SupportFiles\CalibrationData\Spectral Calibration\ParamFiles
2. Select "MtxSTD{Genescan_SetG5}" to open the PAR file.
3. Change "Condition Bounds Range" to [1.0, 20.0].
4. If the calibration was unsuccessful, also change "Sensitivity" to 0.1 and "Quality" to 0.8.
5. Select "Save As" in the "File" menu and save the parameter file under a new name; e.g., MtxStd{Genescan_SetG5_BT5}.par.

Note: Always use this parameter file for spectral calibration runs using QIAGEN Matrix Standard BT5.

Plate Editor for spectral calibration

1. Place the 96-well plate on the autosampler tray.
2. Run the ABI PRISM 3100 Data Collection Software.
3. In Plate View, click "New" to open the Plate Editor dialog box.
4. Enter a name for the plate.
5. Select a Spectral Calibration.
6. Select "96-Well" as plate type, and click "Finish".

Table 16. Plate Editor for spectral calibration

Parameter	Settings
Sample Name	Enter name for the matrix samples
Dye Set	G5
Spectral Run Module	Default (e.g., Spect36_POP4)
Spectral Parameters	MtxStd{GeneScan_SetG5_BT5}.par (parameters created before)

7. Click the column header to select the entire column, and select "Fill Down" from the "Edit" menu to apply the information to the selected samples. Confirm by clicking "OK".

8. Link the reaction plate on the autosampler tray with the created plate ID and start the run.
9. Upon completion of the run, check the Spectral Calibration Result dialog box to show that all capillaries have successfully passed calibration (label A).

If individual capillaries are labeled X, refer to the *ABI PRISM 3100-Avant/3100 Genetic Analyzer User's Manual*.

10. Click "OK" to confirm completion of the run.

Checking the matrix

1. Select "Display Spectral Calibration" from the "Tools" menu, then "Dye Set", and "G5" to review the spectral calibration profile for each capillary.
2. The quality value (Q value) must be greater than 0.95 and the condition number (C value) must be between 1 and 20. Both values must be within the predetermined range.
3. Check for a flat baseline in the matrix samples. There should be 5 peaks with heights of 1000–5000 RFU in each matrix sample.
Note: The optimal range is 2000–4000 RFU.
4. Check the new matrix with the current samples. There should be no pull-up peaks between the dye panels (B, G, Y, R, and O) with the new matrix.
5. If the calibration failed, follow instructions in the section "Spectral parameter" on page 23.
6. If all capillaries have passed the calibration, the last calibration file for Dye Set G5 must be activated manually. Go to "Tools" and click "Set Active Spectral Calibration".
7. Rename the calibration file under "Set Matrix Name" (e.g., BT5_Date of calibration).

Sample preparation

1. Set up a mixture of formamide and DNA size standard according to Table 17.

Table 17. Setup of formamide and DNA size standard mixture

Component	Volume per sample
Hi-Di Formamide	12 μ l
DNA Size Standard 550 (BTO)	0.5 μ l

2. **For each sample to be analyzed, aliquot 12 μ l of the mixture to a tube.**
3. **Add 1 μ l PCR product or allelic ladder (diluted, if necessary).**
4. **Denature for 3 min at 95°C.**
5. **Snap freeze by placing the plate on ice for 3 min.**
Alternatively, the thermal cycler set to 4°C may be used to cool the plate.
6. **Load the samples on the tray.**

Since injections take place simultaneously on all capillaries, 4 or 16 samples must be pipetted onto the plate of multi-capillary analyzers. If fewer samples are analyzed, the empty positions must be filled with 12 μ l Hi-Di Formamide.

To ensure reliable allelic assignment on multi-capillary analyzers, several ladders should be run.

Room temperature may influence the performance of PCR products on multi-capillary instruments, so that shoulder peaks or split peaks occur, especially at low temperatures. Ensure ambient conditions are kept as recommended by the instrument manufacturer.

Setting up the GeneScan Software

1. **Edit the default run module in Dye Set G5 once for the first run.**
Select "Module Editor" to open the dialog box.
2. **Select the appropriate Run Module as template from the GeneScan table (see Table 18 on page 26).**
3. **Modify the Injection Voltage to 3 kV and the Injection Time to 10 s.**
4. **Click "Save As" and enter the name of the new module (e.g., 3kV_10s_260bp). Confirm by clicking "OK".**
5. **Click "Close" to exit the Run Module Editor.**

Table 18. Run Module 3kV_10s_260bp for the ABI PRISM 3100-Avant/3100 Genetic Analyzer

Parameter	Setting
Run Temperature (°C)	Default
Cap Fill Volume	Default
Maximum Current (A)	Default
Current Tolerance (A)	Default
Run Current (A)	Default
Voltage Tolerance (kV)	Default
Pre-Run Voltage (kV)	Default
Pre-Run Time (s)	Default
Injection Voltage (kV)	3.0
Injection Time (s)	10*
Run Voltage (kV)	Default
Number of Steps	Default
Voltage Step Interval	Default
Data Delay Time (s)	Default
Run Time (min)	15 [†]

* Deviating from the standard settings, the injection time may range between 1 and 20 s depending on the type of sample. If samples with very high signal intensities are recorded, a shorter injection time may be selected. For samples with low DNA content, an injection time of up to 20 s may be necessary.

[†] The run time for Investigator DIPlex was modified in order to be able to analyze fragments with lengths of up to 200 bp.

Starting the run

- 1. Place the prepared 96-well plate on the autosampler tray.**
- 2. Run the ABI PRISM 3100 Data Collection Software.**
- 3. In "Plate View", click "New" to open the Plate Editor dialog box.**
- 4. Enter a name for the plate.**
- 5. Select "GeneScan" as the application type.**
- 6. Select "96-Well" as plate type, and click "Finish".**

Table 19. Settings in Plate Editor

Parameter	Settings
Sample Name	Enter name for the matrix samples
Dyes	O
Color Info	Ladder or sample
Project Name	e.g., 3100_Project1
Dye Set	G5
Run Module	3kV_10s_260bp*
Analysis Module 1	DefaultAnalysis.gsp

* See Table 18, "Run Module 3kV_10s_260bp for the ABI PRISM 3100-Avant/3100 Genetic Analyzer".

- 7. Complete the table in the "Plate Editor" and click "OK".**
- 8. Click the column header to highlight the entire column and select "Fill Down" from the "Edit" menu to apply the information to the selected samples.**
- 9. Link the reaction plate on the autosampler tray to the created plate ID and start the run.**
- 10. Upon completion of the run, view the data as Color Data in the Array View of the 3100 Data Collection Software or as Analyzed Sample Files under D:/AppliedBio/3100/DataExtractor/ExtractRuns.**

Analysis parameters

Table 20 lists the recommended analysis parameters.

Table 20. Recommended analysis parameters for the ABI PRISM 3100-Avant/3100 Genetic Analyzer

Parameter	Settings
Analysis Range	Start: 2000 Stop: 6000
Data Processing	Baseline: Checked Multi-component: Checked Smooth options: Light
Peak Detection	Peak Amplitude Thresholds B:* Y:* G:* R:* O:* Min. Peak Half Width: 2 pts Polynomial Degree: 3 Peak Window Size: 11 pts [†]
Size Call Range	Min: 60 Max: 550
Size Calling Method	Local Southern Method
Split Peak Correction	None

* The peak amplitude threshold (cutoff value) corresponds to the minimum peak height that will be detected by the GeneScan or GeneMapper *ID* Software. Thresholds are usually 50–200 RFU and should be determined individually by the laboratory. Recommendation: The minimal peak height should be 3 times higher than the background noise of the baseline.

[†] Only the setting for Peak Window Size is different to defaults from Applied Biosystems for HID analysis.

Note: For information on using the recommended Template Files (as analysis parameters), refer to the appropriate Investigator template files user guide (Genotyper, GeneMapper *ID*, or GeneMapper *ID-X*).

Protocol: Electrophoresis Using the Applied Biosystems 3130/3130xl Genetic Analyzer

For detailed instructions on instrument setup, spectral calibration, or application of the ABI PRISM Data Collection Software version 3.0 and the GeneMapper ID Software, refer to the *Applied Biosystems 3130/3130xl Genetic Analyzers Getting Started Guide*.

The system with 4 capillaries is the Applied Biosystems 3130 Genetic Analyzer, and the system with 16 capillaries is the Applied Biosystems 3130xl Genetic Analyzer.

The virtual filter set Any5Dye is used for combined application of the 5 fluorescent labels 6-FAM, BTG, BTY, BTR, and BTO. This matrix standard is known as BT5.

The materials required for electrophoresis are given in Table 21.

Table 21. Materials needed for electrophoresis

Material	Specifications
Capillary	36 cm Capillary Array for Applied Biosystems 3130/3130xl Genetic Analyzer
Polymer	POP-4 Polymer for Applied Biosystems 3130/3130xl Genetic Analyzer
Buffer	10x Genetic Analyzer Buffer with EDTA

Spectral calibration/matrix generation

Before conducting DNA fragment size analysis, it is necessary to perform a spectral calibration with the 5 fluorescent labels 6-FAM, BTG, BTY, BTR, and BTO for each analyzer. The calibration procedure creates a matrix which is used to correct the overlapping of fluorescence emission spectra of the dyes.

Spectral calibration is comprised of the following steps:

- Preparation of the spectral calibration standards
- Loading the standards to the 96-well reaction plate (one sample per capillary)
- Creating the instrument protocol for spectral calibration (Protocol Manager)
- Defining the plate composition in the plate editor (Plate Manager)
- Performing a spectral calibration run and checking the matrix

Preparing the spectral calibration standards

Example for 4 capillaries (Applied Biosystems 3130 Genetic Analyzer)

1. Set up a mixture of formamide and Matrix Standard BT5 according to Table 22.

Table 22. Setup of formamide and Matrix Standard BT5 mixture for 4 capillaries

Component	Volume
Hi-Di Formamide	60 μ l
Matrix Standard BT5 multi cap.	5 μ l

2. Load 12 μ l of the mixture to a 96-well plate, e.g., positions A1–D1.
3. Denature for 3 min at 95°C.
4. Snap freeze by placing the plate on ice for 3 min.

Alternatively, the thermal cycler set to 4°C may be used to cool the plate.

Example for 16 capillaries (Applied Biosystems 3130xl Genetic Analyzer)

1. Set up a mixture of formamide and Matrix Standard BT5 according to Table 23.

Table 23. Setup of formamide and Matrix Standard BT5 mixture for 16 capillaries

Component	Volume
Hi-Di Formamide	204 μ l
Matrix Standard BT5 multi cap.	17 μ l

2. Load 12 μ l of the mixture to a 96-well plate, e.g., positions A1–H1 and A2–H2.
3. Denature for 3 min at 95°C.
4. Snap freeze by placing the plate on ice for 3 min.

Alternatively, the thermal cycler set to 4°C may be used to cool the plate.

Performing spectral calibration run

1. Place the 96-well plate on the autosampler tray.
2. In the Protocol Manager of the Data Collection Software, open the "Instrument Protocol" window.
3. Click "New" to open the "Protocol Editor" dialog box.
4. Complete the dialog box with information from Table 24 and click "OK".

Table 24. Instrument protocol for spectral calibration

Protocol Editor	Settings
Name	User (e.g., Spectral36_POP4_BT5)
Type	SPECTRAL
Dye Set	Any5Dye
Polymer	User (e.g., POP-4)*
Array Length	User (e.g., 36 cm)*
Chemistry	Matrix Standard
Run Module	Default (e.g., Spect36_POP4_1)*

* Depends on the type of polymer and length of capillary used.

5. Click "New" in the Plate Manager of the Data Collection Software to open the "New Plate" dialog box.
6. Enter information from Table 25 and click "OK". A new table in the Plate Editor opens automatically (Table 26).


Table 25. Plate Editor for spectral calibration (I)

New plate dialog	Settings
Name	e.g., Spectral_BT5_date
Application	Spectral Calibration
Plate Type	96-well
Owner Name/ Operator Name	...

Table 26. Plate Editor for spectral calibration (II)

Parameter	Settings
Sample Name	Enter name for the matrix samples
Priority	e.g., 100
Instrument Protocol 1	Spectral36_POP4_BT5 (setting described before)

- Click the column header to select the entire column, and select “Fill Down” from the “Edit” menu to apply the information to the selected samples. Confirm by clicking “OK”.
- Link the reaction plate on the autosampler tray with the created plate ID (position A or B) and start the run.

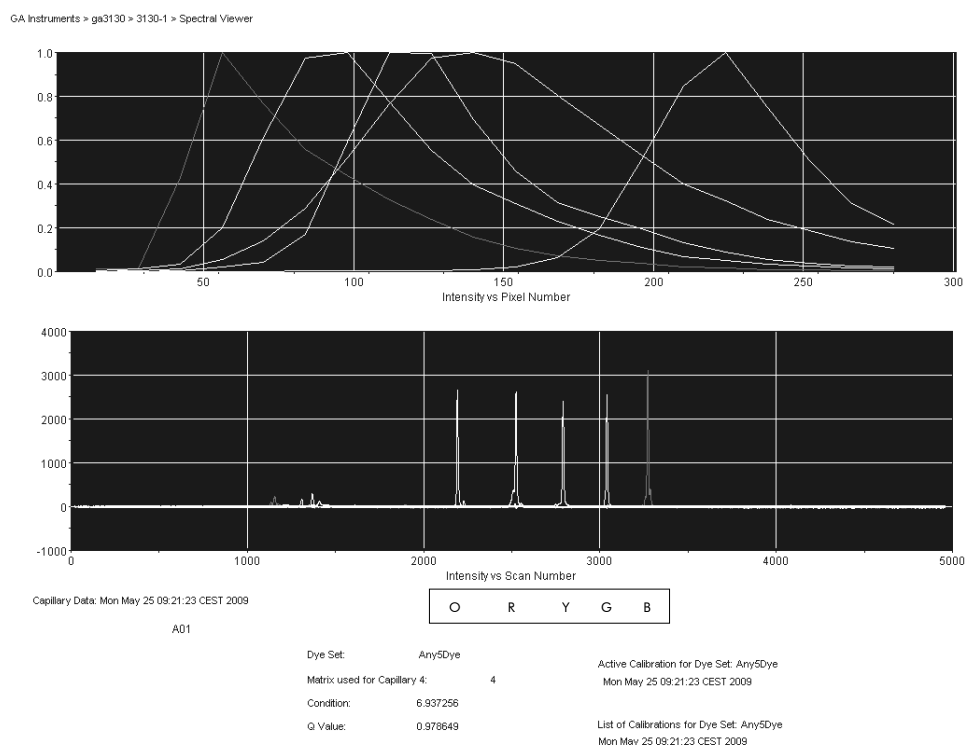

Electropherogram of spectral calibration with matrix standard BT5 on an Applied Biosystems 3130 Genetic Analyzer.

Checking the matrix

1. The quality value (Q value) of each capillary must be greater than 0.95 and the condition number range (C value) must be between 1 and 20.
2. Check for a flat baseline in the matrix samples. As shown in the figure on the previous page, there should be 5 peaks with peak heights of about 1000–5000 RFU in each matrix sample.
Note: The optimal range is 2000–4000 RFU.
3. Check the new matrix with the current samples. There should be no pull-up peaks between the dye panels (B, G, Y, R, O) with the new matrix.
4. If calibration failed, use the optimized values of the Matrix Standard BT5 and repeat the calibration run.
5. If all capillaries have passed the test, the last calibration file for the Dye Set Any5Dye is activated automatically in the Spectral Viewer. Rename the calibration file (e.g., BT5_Date of calibration).

Sample preparation

1. Set up a mixture of formamide and DNA size standard according to Table 27.

Table 27. Setup of formamide and DNA size standard mixture

Component	Volume per sample
Hi-Di Formamide	12.0 μ l
DNA Size Standard 550 (BTO)	0.5 μ l

2. For each sample to be analyzed, aliquot 12 μ l of the mixture to a tube.
3. Add 1 μ l PCR product or allelic ladder (diluted, if necessary).
4. Denature for 3 min at 95°C.
5. Snap freeze by placing the plate on ice for 3 min.
Alternatively, the thermal cycler set to 4°C may be used to cool the plate.
6. Load the samples on the tray.

Since injections take place simultaneously on all capillaries, 4 or 16 samples must be pipetted onto the plate of multi-capillary analyzers. If fewer samples are analyzed, the empty positions must be filled with 12 μ l Hi-Di Formamide.

To ensure reliable allelic assignment on multi-capillary analyzers, several ladders should be run.

Room temperature may influence the performance of PCR products on multi-capillary instruments, so that shoulder peaks or split peaks occur especially at low temperatures. Ensure ambient conditions are kept, as recommended by the instrument manufacturer.

Setting up the Data Collection Software

1. **Edit the Run Module once for the first run. In the Module Manager of the Data Collection Software, click "New" to open the "Run Module Editor" dialog box.**

Note: Modify the Run Module Default settings from "HIDFragmentAnalysis36_POP4_1" to those shown in Table 28

2. **Modify the Injection Voltage to 3 kV and the Injection Time to 10 s (Table 28).**
3. **Click "Save As", enter a name for the new Run Module (e.g., 3kV_10s_260bp), and confirm by clicking "OK".**
4. **Click "Close" to exit the Run Module Editor.**

Table 28. Run Module 3kV_10s_260bp for the Applied Biosystems 3130/3130xl Genetic Analyzer

Parameter	Settings
Oven Temperature (°C)	Default
Poly Fill Volume	Default
Current Stability (μA)	Default
Pre-Run Voltage (kV)	Default
Pre-Run Time (s)	Default
Injection Voltage (kV)	3.0
Injection Time (s)	10*
Voltage Number of Steps	Default
Voltage Step Interval	Default
Data Delay Time (s)	Default
Run Voltage (kV)	Default
Run Time (s)	900 [†]

* Deviating from the standard settings, the injection time may range between 1 and 20 s depending on the type of sample. If samples with very high signal intensities are recorded, a shorter injection time may be selected. For samples with low DNA content, an injection time of up to 20 s may be necessary.

[†] Depending on the analysis conditions the run time for Investigator DIPplex was modified in order to be able to analyze fragments with lengths of up to 200 bp.

Starting the run

- 1. Place the prepared 96-well plate on the autosampler tray.**
- 2. Open the Protocol Manager of the Data Collection Software.**
- 3. Click “New” in the Instrument Protocol window to open the “Protocol Editor” dialog box and enter the information in Table 29.**
- 4. Click “OK” to exit the Protocol Editor.**

Table 29. Settings in Instrument Protocol

Protocol Editor	Settings
Name	Run36_POP4_BT5_15min
Type	REGULAR
Run Module	3kV_10s_260bp*
Dye Set	Any5Dye

* See Table 28, "Run Module 3kV_10s_260bp for the Applied Biosystems 3130/3130xl Genetic Analyzer".

5. Before each run, it is necessary to create a plate definition. In the Plate Manager of the Data Collection Software, click "New" to open the New Plate Dialog box.
6. Enter the information in Table 30.

Table 30. GeneMapper Plate Editor (I)

Protocol Editor	Settings
Name	e.g., Plate_BT5_Date
Application	Select GeneMapper Application
Plate type	96-Well
Owner Name/ Operator Name	...

7. Click "OK" and a new table in the Plate Editor opens automatically (Table 31).
8. Click the column header to select the entire column. Select "Fill Down" from the "Edit" menu to apply the information to all selected samples. Click "OK".
9. In the Run Scheduler, click "Find All", and select "Link" to link the reaction plate on the autosampler tray to the newly created plate record (position A or B).

Table 31. GeneMapper Plate Editor (II)

Parameter	Settings
Sample Name	Enter the name for the samples
Priority	e.g., 100 (Default)
Sample Type	Sample or Allelic Ladder
Size Standard	e.g., SST-BTO_60-260bp
Panel	e.g., DIPplex_Panels
Analysis Method	e.g., DIPplex_HID_3130_200rfu
Snp Set	–
User-defined 1-3	–
Results Group 1	(Select results group)
Instrument Protocol 1	Run36_POP4_BT5_15min (setting described before)

10. Start the run.**11. During the run, view Error Status in the Event Log or examine the quality of the raw data for each capillary in the Capillaries Viewer or the Cap/Array Viewer.****12. View data as an overview in Run History or Cap/Array Viewer of the Data Collection Software.**

Run data are saved in the Run Folder of the previously chosen Result Group.

Analysis parameters/analysis method

Table 32 lists the recommended analysis parameters in the worksheet “Peak Detector”. Table 33 and Table 34 list the recommended settings for the worksheet “Allele” and “Peak Quality”.

Table 32. Recommended settings for the Applied Biosystems 3130/3130xl Genetic Analyzer

Parameter	Settings
Peak Detection Algorithm	Advanced
Ranges	Analysis: Partial Range Start Point: 2000; Stop Point: 6000 Sizing: All Sizes
Smoothing and Baseline	Smoothing: Light Baseline Window: 51 pts
Size Calling Method	Local Southern Method
Peak Detection	Peak Amplitude Thresholds B:* Y:* G:* R:* O:* Min. Peak Half Width: 2 pts Polynomial Degree: 3 Peak Window Size: 11 pts [†] Slope Thresholds: 0.0

* The peak amplitude threshold (cutoff value) corresponds to the minimum peak height that will be detected by the GeneMapper *ID* Software. The thresholds are usually 50–200 RFU and should be determined individually by the laboratory. Recommendation: The minimal peak height should be 3 times higher than the background noise of the baseline.

[†] Only the setting for Peak Window Size is different to defaults from Applied Biosystems for HID analysis.

Table 33. Recommended settings in the worksheet “Allele”

Parameter	Settings
Amelogenin Cutoff	0.1*

* As with Amelogenin, all DIPs will be examined using GeneMapper *ID* Software.

Table 34. Recommended settings in the worksheet “Peak Quality”

Parameter	Settings
Heterozygote Balance (Minimum peak height ratio)	0.1
Allele Number (Maximum expected alleles)	20

Note: For information on using the recommended Template Files (as analysis parameters), refer to the appropriate Investigator Template files user guide (Genotyper, GeneMapper *ID*, or GeneMapper *ID-X*).

Protocol: Electrophoresis Using the Applied Biosystems 3500/3500xL Genetic Analyzer

For detailed instructions on instrument setup, spectral calibration, or application of the Applied Biosystems 3500 Series Data Collection Software version 1.0 and the GeneMapper ID-X Software version 1.2, refer to the *Applied Biosystems 3500/3500xL Genetic Analyzers User Guide*.

The system with 8 capillaries is the Applied Biosystems 3500 Genetic Analyzer and the system with 24 capillaries is the Applied Biosystems 3500xL Genetic Analyzer.

The virtual filter set AnyDye is used for combined application of the 5 fluorescent labels 6-FAM, BTG, BTY, BTR, and BTO. This matrix standard is known as BT5.

The materials required for electrophoresis are given in Table 35.

Table 35. Materials required for electrophoresis

Material	Specifications
Capillary	36 cm Array for Applied Biosystems 3500/3500xL Genetic Analyzer
Polymer	POP-4 for Applied Biosystems 3500/3500xL Genetic Analyzer
Buffer	Anode Buffer Container (ABC) 3500 Series Cathode Buffer Container (CBC) 3500 Series

Spectral calibration/matrix generation

Before conducting DNA fragment size analysis, it is necessary to perform a spectral calibration with the 5 fluorescent labels 6-FAM, BTG, BTY, BTR, and BTO for each analyzer. The calibration procedure creates a matrix which is used to correct the overlapping of fluorescence emission spectra of the dyes.

IMPORTANT: Spectral calibration must be performed for each new capillary array.

Spectral calibration is comprised of the following steps:

- Preparation of the instrument
- Preparation of dye set BT5
- Preparation of the standard calibration plate
- Plate assembly and loading the plate in the instrument

- Performing a spectral calibration run
- Checking the matrix

Preparation of the instrument

Before the spectral calibration process, ensure that the spatial calibration has been performed. This process is described in detail in the *Applied Biosystems 3500/3500xL Genetic Analyzers User Guide*.

Table 36. The fluorescent labels of BT5

Color	Matrix standard
Blue (B)	6-FAM
Green (G)	BTG
Yellow (Y)	BTY
Red (R)	BTR
Orange (O)	BTO

Software setup of dye set BT5

Before the spectral calibration, a dye set for the Matrix Standard BT5 must be set up.

1. To create a new dye set, go to “Library” and select “Analyze”, followed by “Dye Sets”, and click “Create”.
2. Enter a “Dye Set Name”, e.g., BT5.
3. Select “Matrix Standard” as a chemistry and “AnyDye Template” as a dye set template.
4. Disable “Purple” in the field “Arrange Dyes”. Ensure that all other colors are enabled.
5. Under “Calibration Peak Order” the colors need to be arranged as follows: 5 — blue, 4 — green, 3 — yellow, 2 — red, and 1 — orange.
6. Do not change the “Parameter” settings.
7. Click “Save” to confirm the changes.

Create New Dye Set

Setup a Dye Set

* Dye Set Name: BT5

* Chemistry: Matrix Standard

* Dye Set Template: AnyDye Template

Arrange Dyes

Dye Selection	Reduced Selection	Calibration Peak Order
☒	☒	5
☒	☒	4
☒	☒	3
☒	☒	2
☐	☒	0
☒	☒	1

Parameters

The parameters will be used for instruments configured with 36cm capillary array and polymer POP4

Matrix Condition Number Upper Limit: 20.0

Locate Start Point: * After Scan 300 * Before Scan 5000

* Limit Scans To: 20000

Sensitivity: 0.1

* Minimum Quality Score: 0.8

Notes

Matrix Std. BT5 multi cap.

Close Save

Setup of dye set BT5.

Preparation of the standard calibration plate

Example for 8 capillaries (Applied Biosystems 3500 Genetic Analyzer)

1. Set up a mixture of formamide and Matrix Standard BT5 according to Table 37.

Table 37. Setup of formamide and Matrix Standard BT5 mixture for 8 capillaries

Component	Volume
Hi-Di Formamide	90 μ l
Matrix Standard BT5 multi cap.	10 μ l

2. Load 10 μ l of the mixture to a 96-well plate, e.g., positions A1–H1.
3. Denature for 3 min at 95°C.
4. Snap freeze by placing the plate on ice for 3 min.

Alternatively, the thermal cycler set to 4°C may be used to cool the plate.

Example for 24 capillaries (Applied Biosystems 3500xL Genetic Analyzer)

1. Set up a mixture of formamide and Matrix Standard BT5 according to Table 38.

Table 38. Setup of formamide and Matrix Standard BT5 mixture for 24 capillaries

Component	Volume
Hi-Di Formamide	225 μ l
Matrix Standard BT5 multi cap.	25 μ l

2. Load 10 μ l of the mixture to a 96-well plate, e.g., positions A1–H1, A2–H2, and A3–H3. When using a 384-well plate, load 10 μ l of the mixtures to columns 1, 3, and 5 in rows A, C, E, G, I, K, M, and O.
3. Denature for 3 min at 95°C.
4. Snap freeze by placing the plate on ice for 3 min.

Alternatively, the thermal cycler set to 4°C may be used to cool the plate.

Plate assembly and loading the plate in the instrument

The necessary steps are described in detail in the *Applied Biosystems 3500/3500xL Genetic Analyzers User Guide*.

Performing a spectral calibration run


Once the multiwell plates containing the spectral calibration mixture are placed in the autosampler tray, the spectral calibration process can be started.

1. To access the Spectral Calibration screen, select “Maintenance” on the Dashboard of the 3500 Series Data Collection software.
2. The number of wells in the spectral calibration plate and their location in the instrument must be specified.
3. Select “Matrix Standard” as a chemistry standard and “BT5” for dye set.
(Optional) Enable “Allow Borrowing”.
4. Click “Start Run”.

Checking the matrix

1. Click a capillary in the table to display the results for each capillary (spectral data, quality value, and condition number) below the run results table.
2. The quality value (Q value) of each capillary must be greater than 0.8 and the number range (C value) must be between 1 and 20.
3. Check the matrix samples for a flat baseline. As shown in the figure, there should be 5 peaks with peak heights of about 1000–5000 RFU for each matrix sample

Note: The optimal range is 2000–4000 RFU.

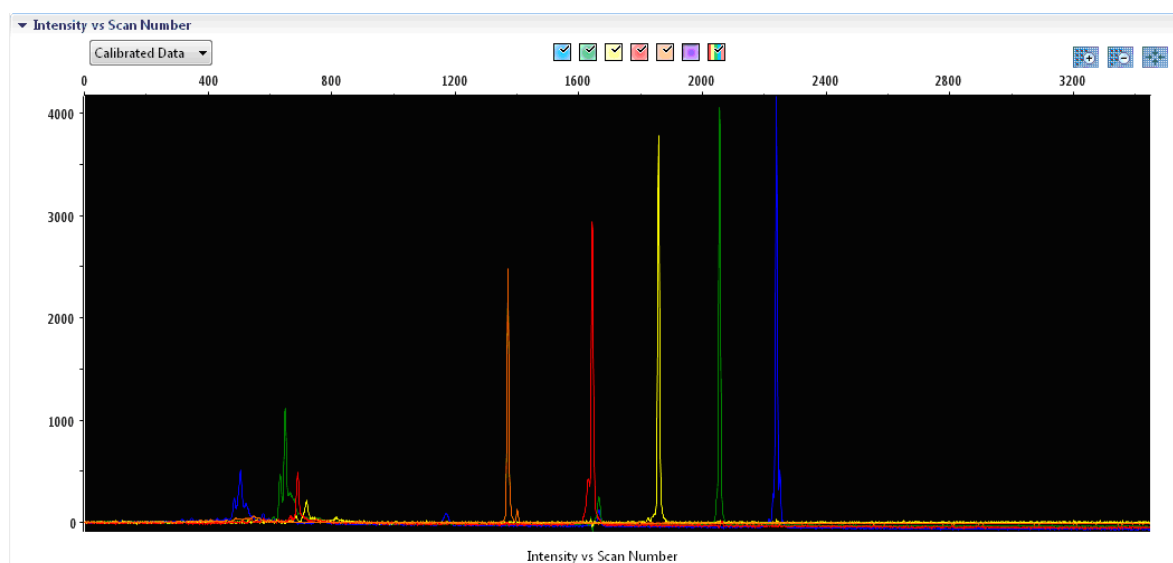

Electropherogram of spectral calibration with matrix standard BT5 on an Applied Biosystems 3500 Genetic Analyzer.

4. If all capillaries have passed the calibration, the “Overall” row displays green results. If the “Overall” row displays red results, refer to the “Spectral Calibration Troubleshooting” section of the *Applied Biosystems 3500/3500xL Genetic Analyzers User Guide*.

▼ Capillary Run Data								
Capillary	1	2	3	4	5	6	7	8
Run 1	Passed	Passed	Passed	Passed	Passed	Passed	Passed	Passed
Run 2	Not Calibrated	Not Calibrated	Not Calibrated	Not Calibrated	Not Calibrated	Not Calibrated	Not Calibrated	Not Calibrated
Run 3	Not Calibrated	Not Calibrated	Not Calibrated	Not Calibrated	Not Calibrated	Not Calibrated	Not Calibrated	Not Calibrated
Overall	Passed	Passed	Passed	Passed	Passed	Passed	Passed	Passed

■ Passed
 ■ Failed
 ■ Borrowed
 ■ Not Calibrated

Example of a successful spectral calibration of the matrix standard BT5 for all capillaries on an Applied Biosystems 3500 Genetic Analyzer.

5. **For each capillary, select and display the spectral and raw data. Check that the data meet the following criteria:**

The order of the peaks in the spectral profile from left to right read orange-red-yellow-green-blue.

No extraneous peaks appear in the raw data profile.

Peak morphology in the spectral profile shows no gross overlaps, dips, or other irregularities. Separate and distinct peaks should be visible.

6. **If the data for all capillaries meet the criteria above, click "Accept Results". If any capillary data does not meet the criteria above, click "Reject Results", and refer to the "spectral calibration troubleshooting" section of the *Applied Biosystems 3500/3500xL Genetic Analyzers User Guide*.**

Sample preparation

1. **Set up a mixture of formamide and DNA size standard according to Table 39.**

Table 39. Setup of formamide and DNA size standard mixture

Component	Volume per sample
Hi-Di Formamide	12.0 μ l
DNA Size Standard 550 (BTO)	0.5 μ l

2. **For each sample to be analyzed, aliquot 12 μ l of the mixture to a tube.**
3. **Add 1 μ l PCR product or allelic ladder (diluted, if necessary).**
4. **Denature for 3 min at 95°C.**
5. **Snap freeze by placing the plate on ice for 3 min.**
Alternatively, the thermal cycler set to 4°C may be used to cool the plate.
6. **Load the samples on the tray.**

Since injections take place simultaneously on all capillaries, 8 or 24 samples must be pipetted onto the plate of multi-capillary analyzers. If fewer samples are analyzed, the empty positions must be filled with 12 μ l Hi-Di Formamide.

To ensure a reliable allelic assignment on multi-capillary analyzers, inject one allelic ladder for each set of 24 samples:

- 8-capillary instruments: One allelic ladder per 3 injections
- 24-capillary instruments: One allelic ladder per 1 injection

Room temperature may influence the performance of PCR products on multi-capillary instruments, so that shoulder peaks or split peaks occur, especially at low temperatures. Ensure ambient conditions are kept, as recommended by the instrument manufacturer.

Setting up a run

If using the Investigator DIPplex Kit for the first time on an Applied Biosystems 3500 Genetic Analyzer, you will first need to set up a number of protocols:

- Instrument Protocol
- Size Standard
- QC Protocol
- Assay

All protocols can be set up via the Dashboard of the 3500 Series Data Collection software.

- 1. To set up the Instrument Protocol, go to “Library” and select “Analyze”, followed by “Instrument Protocols”, and click “Create”.**

Note: Modify the Run Module Default settings from “HID36_POP4” as shown in Table 40.

- 2. The parameters in Table 40 must be entered or selected.**

Table 40. Instrument Protocol parameters for Applied Biosystems 3500/3500xL Genetic Analyzer

Parameter	Setting
Application Type	HID
Capillary Length	36 cm
Polymer	POP-4
Dye Set	e.g., BT5
Run Module	HID36_POP4
Protocol Name	e.g., Investigator DIPplex
Oven Temperature (°C)	Default
Run Voltage (kV)	Default
Pre-Run Voltage (kV)	Default
Injection Voltage (kV)	3.0
Run Time (s)	900
Pre-Run Time (s)	Default
Injection Time (s)	8.0*
Data Delay (s)	Default
Advanced Options	Default

* Deviating from the standard settings, the injection time may range between 1 and 20 s depending on the type of sample. If samples with very high signal intensities are recorded, a shorter injection time may be selected. For samples with low DNA content, an injection time of up to 20 s may be necessary.

3. Click **"Save"** to confirm the changes.
4. To set up the Size Standard, go to **"Library"**, select **"Analyze"**, followed by **"Size Standards"**, and click **"Create"**.
5. The parameters in Table 41 must be entered or selected.

The DNA Size Standard 550 (BTO) should be used with the following lengths of fragments: 60, 80, 90, 100, 120, 140, 160, 180, 200, 220, 240, 250, and 260 bp.

Table 41. Size standard parameters

Parameter	Setting
Size Standard	e.g., SST-BTO_60-260bp
Dye Color	Orange

6. Click **"Save"** to confirm the changes.
7. To set up the QC Protocol, go to **"Library"** and select **"Analyze"**, followed by **"QC Protocols"**, and click **"Create"**.
8. The parameters in Table 42 must be entered or selected.

Table 42. QC Protocol parameters

Parameter	Setting
Protocol Name	e.g., BTO_550
Size Standard	SST-BTO_60-260 (from Table 41)
Sizecaller	SizeCaller v1.1.0

9. Go to **"Analysis Settings"**, followed by **"Peak Amplitude Threshold"** and disable **"Purple"**. Ensure that all other colors are enabled.
All other settings should remain as **"Default"**.
10. Click **"Save"** to confirm the changes.
11. To set up an Assay, go to **"Library"** and select **"Manage"**, followed by **"Assays"**, and click **"Create"**.
12. To analyze Investigator DIPplex fragments, the parameters in Table 43 must be selected.

Table 43. Assay parameters

Parameter	Setting
Assay Name	e.g., Investigator DIPplex
Color	Default
Application Type	HID
Instrument Protocol	e.g., Investigator DIPplex (from Table 40)
QC Protocols	e.g., BTO_550 (from Table 42)

13. Click “Save” to confirm the changes.

Starting the run

- 1. In the Dashboard, click “Create New Plate”.**
- 2. Go to “Define Plate Properties” and select “Plate Details”. Select or enter the parameters in Table 44.**

Table 44. Plate properties

Property	Setting
Name	e.g., Investigator DIPplex
Number of Wells	96
Plate Type	HID
Capillary Length	36 cm
Polymer	POP-4

- 3. Click “Assign Plate Contents” to confirm the changes.**
- 4. Enter the designated sample name in each well containing a sample or allelic ladder. This will identify the well positions of each sample for the data collection and processing.**
- 5. Choose the correct Assay for the analysis. If you followed the steps under “Setting up a run” (page 46), this would be Investigator DIPplex from step 11. All named wells on the plate must have an assigned assay.**
- 6. Select the wells for which to specify an assay. Check the box next to the assay name to assign it to the selected wells.**

7. (Optional) Repeat for file name conventions and results group.
8. If not already done, load the instrument with the assembled plate and close the instrument door to re-initialize the instrument. Then click "Link Plate for Run". In the next screen, enter the desired Run Name and click "Start Run".

Analysis parameters/analysis method

Table 45 lists the recommended analysis parameters in the worksheet "Peak Detector". Table 46 and Table 47 list the settings for the worksheet "Allele" and "Peak Quality".

Table 45. Recommended settings for the Applied Biosystems 3500/3500xL Genetic Analyzer

Parameter	Settings
Peak Detection Algorithm	Advanced
Ranges	Analysis: Partial Range Start Point: 1000; Stop Point: 20,000 Sizing: All Sizes
Smoothing and Baselineing	Smoothing: Light Baseline Window: 51 pts
Size Calling Method	Local Southern Method
Peak Detection	Peak Amplitude Thresholds B:* Y:* G:* R:* O:* Min. Peak Half Width: 2 pts Polynomial Degree: 3 Peak Window Size: 11 pts [†] Slope Thresholds: 0.0

* The peak amplitude threshold (cutoff value) corresponds to the minimum peak height that will be detected by the GeneMapper *ID-X* Software version 1.2. The thresholds are usually 50–200 RFU and should be determined individually by the laboratory. Recommendation: The minimal peak height should be 3 times higher than the background noise of the baseline.

[†] Only the setting for Peak Window Size is different to defaults from Applied Biosystems for HID analysis.

Table 46. Recommended settings in the worksheet “Allele”

Parameter	Settings
Amelogenin Cutoff	0.1*

* As with Amelogenin, all DIPs will be examined using GeneMapper ID-X Software.

Table 47. Recommended settings in the worksheet “Peak Quality”

Parameter	Settings
Heterozygote Balance (Minimum peak height ratio)	0.1
Allele Number (Maximum expected alleles)	20

Note: For information on using the recommended Template Files (as analysis parameters), refer to the appropriate Investigator template files user guide (Genotyper, GeneMapper ID, or GeneMapper ID-X).

Protocol: Analysis

For general instructions on automatic sample analysis, refer to the *GeneScan*, *GeneMapper ID* or *GeneMapper ID-X Software User Guide*.

Finding the exact lengths of the amplified products depends on the device type, the conditions of electrophoresis, as well as the DNA size standard used. Due to the complexity of some loci, determining the size should be based on evenly distributed references. The DNA Size Standard 550 (BTO) should be used with the following lengths of fragments: 60, 80, 90, 100, 120, 140, 160, 180, 200, 220, 240, 250, 260, 280, 300, 320, 340, 360, 380, 400, 425, 450, 475, 500, 525, and 550 bp.

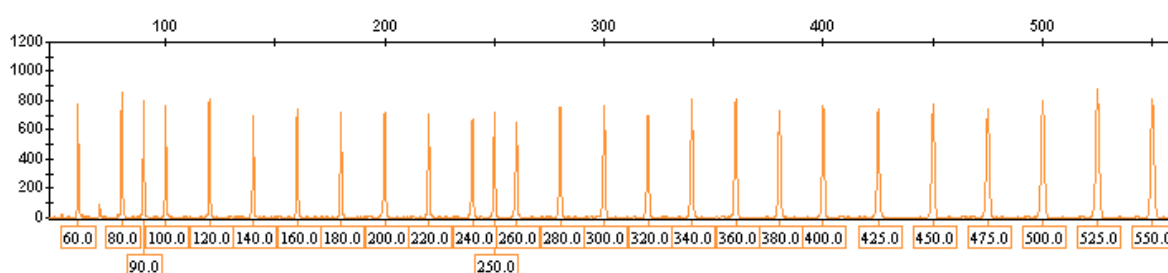


Figure 1. Electropherogram of the DNA Size Standard 550 (BTO), fragments with lengths in bp.

Analysis software

Allele allocation should be carried out with suitable analysis software, e.g., *GeneMapper ID-X*, *GeneMapper ID*, or *Genotyper Software* in combination with the Investigator Template Files available as a download from www.qiagen.com, see Table 48 and Table 49.

The recommended Investigator Template File for *Genotyper Software* is the DIPlex.

In order to sort genotypes for easy interpretation of the results, QIAGEN offers DIPSorter Freeware. If using DIPSorter Freeware, use the recommended Table Settings for displaying and exporting the genotype table from *GeneMapper ID* (Table for 20 Alleles) or *Genotyper software* (Vertical Table for 20 Alleles).

Table 48. Recommended Investigator Template Files for GeneMapper ID

File type	File name
Panels	DIPplex_Panels
BinSets	DIPplex_Bins
Size standard	SST-BTO_60–260bp
Analysis Method	DIPplex_HID_310_50rfu DIPplex_HID_3130_50rfu DIPplex_HID_310_200rfu DIPplex_HID_3130_200rfu
Plot Settings	Plots_5dyes
Table Settings	Table for 20 alleles

Panels and BinSets must always be used; other template files are optional.

Table 49. Recommended Investigator Template Files for GeneMapper ID-X

File type	File name
Panels	DIPplex_Panels_x
BinSets	DIPplex_Bins_x
Size standard	SST-BTO_60–260bp
Analysis Method	DIPplex_HID_310_50rfu DIPplex_HID_3130_50rfu DIPplex_HID_310_200rfu DIPplex_HID_3130_200rfu DIPplex_HID_3500
Plot Settings	Plots_5dyes
Table Settings	Table for 20 alleles

Panels and BinSets must always be used; other template files are optional.

Controls

The alleles listed in Table 50 represent the Control DNA 9948 (included in the Investigator DIPplex Kit) and DNA from other commercially available standard cell lines.

Table 50. Allele assignment of the Investigator DIPplex Kit

Locus	Control DNA 9948	ATCC K-562	CCR 9947A	CCR 3657
Amelogenin	X/Y	X	X	X/Y
HLD77	+	-/+	-	+
HLD45	+	-	-/+	-/+
HLD131	-/+	-/+	-	+
HLD70	-/+	-/+	+	-
HLD6	-/+	+	+	-/+
HLD111	-/+	+	-/+	-/+
HLD58	+	-	+	+
HLD56	-	-	-/+	-/+
HLD118	-	-	-	-/+
HLD92	+	+	-/+	-
HLD93	-	+	-	-/+
HLD99	+	+	-/+	-/+
HLD88	-/+	-	-	+
HLD101	-/+	-/+	-/+	-/+
HLD67	+	-/+	+	+
HLD83	-	-	-/+	-
HLD114	+	-	-	-/+
HLD48	-/+	+	+	+
HLD124	+	-	-/+	+

Table continued on next page.

Table continued from previous page.

Locus	Control DNA 9948	ATCC K-562	CCR 9947A	CCR 3657
HLD122	–	–	–/+	–/+
HLD125	+	–	–/+	–
HLD64	–/+	–	+	–/+
HLD81	–/+	–/+	–	+
HLD136	+	+	–/+	–/+
HLD133	+	–	+	–/+
HLD97	–/+	–	–/+	+
HLD40	–/+	+	–	+
HLD128	–	–/+	–/+	–/+
HLD39	+	–/+	–	–
HLD84	–/+	+	–	–

For further confirmation, the table above displays the alleles of the reference DNA purchased from Coriell Cell Repositories (CCR), as well as 3 reference DNAs purchased from CCR and ATCC up to the standard of Szibor et al. (2).

Lengths of fragments and alleles

Table 51 and Table 52 show the fragment lengths of individual alleles that refer to the DNA Size Standard 550 (BTO). Analyses were performed using an Applied Biosystems 3500 Genetic Analyzer (Figure 2) and an Applied Biosystems 3130 Genetic Analyzer (Figure 3), both with POP-4 polymer. Different analysis instruments, DNA size standards, or polymers may result in different fragment lengths. In addition, a visual alignment with the allelic ladder is recommended.

Scaling

- Horizontal: 70–165 bp
- Vertical: Depending on signal intensity

Table 51. Fragment lengths of the allelic ladder DIPplex analyzed on an ABI PRISM 3130 Genetic Analyzer (blue and green panels)

Marker (blue)	–DIP (bp)*	+DIP (bp)*	Marker (green)	–DIP (bp)*	+DIP (bp)*
Amelogenin	77 (X)	80 (Y)	HLD118	77	81
HLD77	84	88	HLD92	87	90
HLD45	92	96	HLD93	98	103
HLD131	100	113	HLD99	108	113
HLD70	104	108	HLD88	118	128
HLD6	118	134	HLD101	131	135
HLD111	122	126	HLD67	140	147
HLD58	136	140			
HLD56	144	149			

* Rounded to integer.

Table 52. Fragment lengths of allelic ladder DIPplex analyzed on an ABI PRISM 3130 Genetic Analyzer (yellow and red panels)

Marker (yellow)	–DIP (bp)*	+DIP (bp)*	Marker (red)	–DIP (bp)*	+DIP (bp)*
HLD83	76	80	HLD136	79	84
HLD114	83	100	HLD133	91	102
HLD48	89	94	HLD97	97	110
HLD124	104	109	HLD40	106	129
HLD122	115	126	HLD128	113	121
HLD125	129	135	HLD39	125	143
HLD64	138	144	HLD84	135	140
HLD81	147	158			

* Rounded to integer.

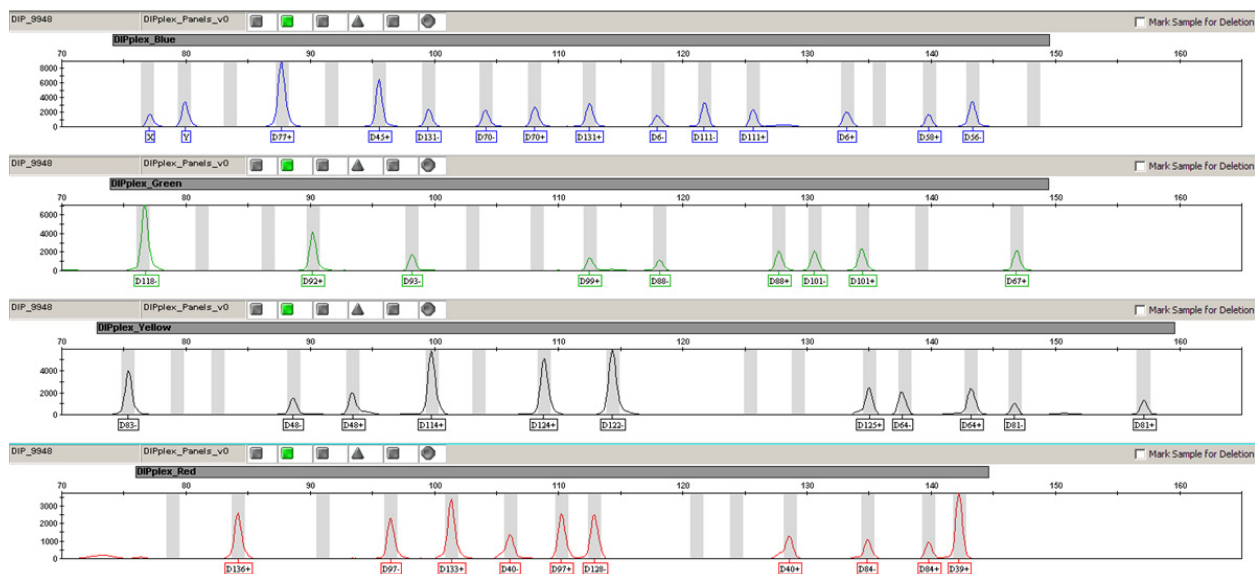


Figure 2. Electropherogram of the Investigator DIPplex Kit using 500 pg Control DNA 9948. Analysis was performed on an Applied Biosystems 3500 Genetic Analyzer. Allele assignment was performed using GeneMapper ID-X Software.

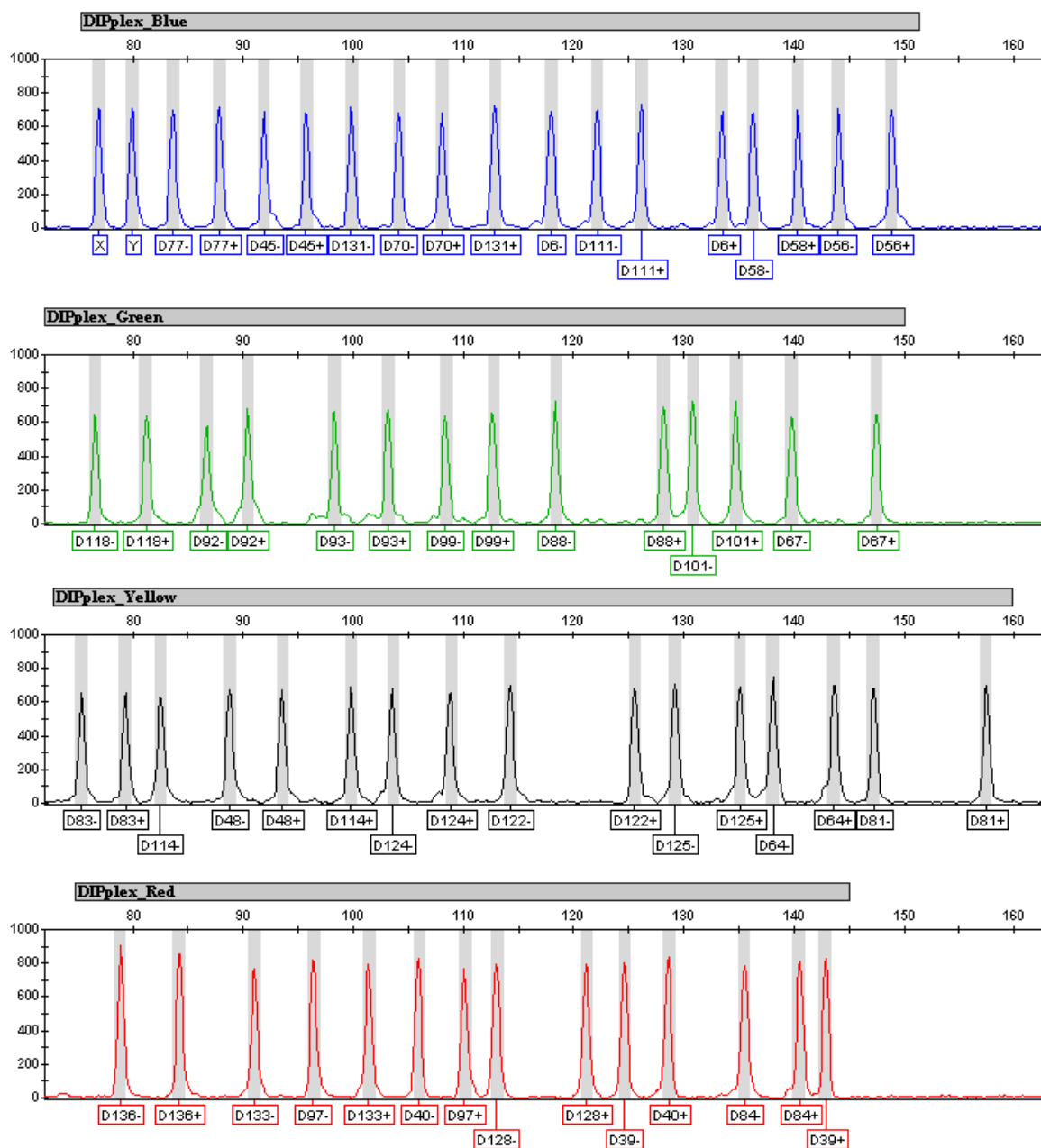


Figure 3. Electropherogram of the allelic ladder DIPplex analyzed on an Applied Biosystems 3130 Genetic Analyzer. Allele assignment was performed using GeneMapper ID Software and the Investigator DIPplex Template File.

Interpretation of Results

Post-PCR analysis and automatic allele assignment with suitable analysis software ensure a precise and reliable discrimination of alleles.

General procedure for the analysis

1. Check the DNA size standard.
2. Check the allelic ladder.
3. Check the positive control.
4. Check the negative control.
5. Analyze and interpret the sample data.

Pull-up peaks

Pull-up peaks may occur if peak heights are outside the linear detection range (see "Troubleshooting Guide", page 60), or if an incorrect matrix was applied. They appear at positions of specific peaks in other color channels, typically with lower signal intensities. Peak heights should not exceed thresholds in order to prevent pull-up peaks.

Template-independent addition of nucleotides

Because of its terminal transferase activity, the *Taq* DNA polymerase may cause incomplete adenylation at the 3'-end of the amplified DNA fragments. The artifact peak is one base shorter than expected (–1 peaks). All primers included in the Investigator DIPplex Kit are designed to minimize these artifacts. Artifact formation is further reduced by the final extension step of the PCR protocol at 68°C for 60 minutes. Peak height of the artifact correlates with the amount of DNA. Laboratories should define their own limits for analysis of the peaks.

Artifacts

Room temperature may influence the performance of PCR products on multi-capillary instruments, so that shoulder peaks or split peaks occur. If shoulder or split peaks appear, we recommend injecting the sample again.

Outside of the allelic range (> 170 bp) amplification artifacts with very low signal intensity may occur.

Troubleshooting Guide

This troubleshooting guide may be helpful in solving any problems that may arise. For more information, see also the Frequently Asked Questions page at our Technical Support Center: www.qiagen.com/FAQ/FAQList.aspx. The scientists in QIAGEN Technical Services are always happy to answer any questions you may have about either the information and protocols in this handbook or sample and assay technologies (for contact information, see back cover or visit www.qiagen.com).

Comments and suggestions

Sample preparation

Sample signal intensity must be increased	Reduce the volume of the DNA Size Standard 550 (BTO) to peak heights of about 500 RFU. Purify the PCR products before starting the analysis. We recommend the MinElute® PCR Purification Kit for rapid and effective purification (see Ordering Information).
---	--

Matrix/spectral calibration is not appropriate

There are pull-up peaks between the dye panels (B, G, Y, R, O) with the current matrix/spectral calibration	This matrix cannot be used for the analysis. Repeat the matrix generation/spectral calibration. Be sure to carefully follow the correct protocol for the specific analysis instrument.
---	--

Many peaks are labeled as off-ladder (OL) alleles in the samples

a) DNA Size Standard 550 (BTO) was not defined or identified correctly	Click the orange "Size Match Editor" icon in the upper toolbar or the GeneMapper ID or GeneMapper ID-X Software. Mark the orange fragments of all samples. Always use the DNA Size Standard 550 included in Investigator Human Identification PCR Kits.
--	--

Comments and suggestions

- | | |
|---|--|
| b) Signal intensities are too high. If the peak heights of the samples are outside the linear detection range (>4000 RFU/>5000 RFU*), stutters, split peaks, and artifacts may be increased | Reduce the injection time in increments to a minimum of 1 s, reduce the amount of the PCR amplification product for analysis, or reduce the quantity of DNA for PCR. |
| c) Bubbles in the capillary lead to pull-up peaks in all color panels ("spikes") that result in allele misnomer | Repeat electrophoresis to confirm results. |
| d) Differences in the run performance among the capillaries of a multi-capillary analyzer may result in allelic assignment shift | For reliable allelic assignment on multi-capillary analyzers, a number of allelic ladders should be run. |
| e) Use of 32-cycle PCR program for small amounts of DNA | Too small amounts of DNA may result in allelic dropouts and imbalances of the peaks. Furthermore, nonspecific amplification products may appear. By increasing the number of cycles, there is a risk of cross-contamination due to impurities. |

Injection/file of the allelic ladder is not appropriate

- | | |
|--|---|
| a) An additional signal can be identified as peak of the allelic ladder because of dysfunctions during the electrophoresis. If peaks of the allelic ladder are miscalled, the ladder cannot be used for the analysis | Use a different injection/file of the allelic ladder and check the data of the analyzed sizes from the Size Standard (in bp) of the allelic ladder.

Always use the DNA Size Standard 550 for Investigator Human Identification PCR Kits. |
|--|---|

* >4000 RFU for the ABI PRISM 310 Genetic Analyzer; >5000 RFU for the ABI PRISM 3100 and Applied Biosystems 3130/3500 Genetic Analyzers.

Comments and suggestions

- | | |
|---|---|
| b) One peak of the allelic ladder is below the peak detection value (50–200 RFU) of the analysis method used, and thus, is not identified | The allelic ladder must be loaded onto the analysis instrument at a higher concentration than samples to be analyzed.

Alternatively, allelic ladder data can be analyzed with a lower peak detection value in Analysis Software. |
| c) One peak of the allelic ladder is not identified because it is outside the expected size range of the software (in bp) | Compare the length of the fragments (in bp) of the first allele in one color of the allelic ladder with the corresponding value in the categories. Then compare it with the other alleles. |
| d) Point alleles are not found | Point alleles are alleles with at least 1 bp difference to the next integer allele. Check the settings of the analysis method. Lower the Peak Window Size value to 11 points. |

References

QIAGEN maintains a large, up-to-date online database of scientific publications utilizing QIAGEN products. Comprehensive search options allow you to find the articles you need, either by a simple keyword search or by specifying the application, research area, title, etc.

For a complete list of references, visit the QIAGEN Reference Database online at www.qiagen.com/RefDB/search.asp or contact QIAGEN Technical Services or your local distributor.

Cited references

1. Sanchez, J.J., et al. (2006) A multiplex assay with 52 single nucleotide polymorphisms for human identification. *Electrophoresis* **27**, 1713.
2. Szibor, R., Krawczak, M., Hering, S., Edelmann, J., Kuhlisch, E., and Krause, D. (2003) Use of X-linked markers for forensic purposes. *Int. J. Legal Med.* **117**, 67.

General references

Bär, W., et al. (1997) DNA recommendations. Further report of the DNA Commission of the ISFG regarding the use of short tandem repeat systems. *Int. J. Legal Med.* **110**, 175.

Gill, P., et al. (2006) DNA Commission of the International Society of Forensic Genetics (ISFG): Recommendations on the interpretation of mixtures. *Forensic Sci. Int.* **160**, 90.

Mills, R.E., et al. (2006) An initial map of insertion and deletion (INDEL) variation in the human genome. *Genome Res.* **16**, 1182.

Weber, J.L., David, D., Heil, J., Fan, Y., Zhao, C., and Marth, G. (2002) Human diallelic insertion/deletion polymorphisms. *Am. J. Hum. Genet.* **71**, 854.

Ordering Information

Product	Contents	Cat. no.
Investigator DIPplex Kit (100)	Primer mix, reaction mix, DNA Polymerase, Control DNA, allelic ladder, DNA size standard, and nuclease-free water	384015
Related products		
Investigator Quantification Kits		
Investigator Quantiplex Kit (200)	Primer mix IC FQ, reaction mix FQ, control DNA Z1, dilution buffer	387016
Investigator Quantiplex HYres Kit	Primer mix IC YQ, reaction mix YQ, control DNA Z1, dilution buffer	387116
Accessories		
DNA Size Standard 550 (BTO) (100)	DNA Size Standard 550 (BTO) for 100 reactions	386015
Matrix Standard BT5 single cap. (5 x 25)	Matrix standard 6-FAM, BTG, BTY, BTR, and BTO	386113
Matrix Standard BT5 multi cap. (25)	Matrix standard 6-FAM, BTG, BTY, BTR, and BTO	386123
Matrix Standard BT5 multi cap. (50)	Matrix standard 6-FAM, BTG, BTY, BTR, and BTO	386125
Multi Taq2 DNA Polymerase (100)	100 Units Multi Taq2 DNA Polymerase	386315

Product	Contents	Cat. no.
DNA extraction and purification		
QIAamp® DNA Investigator Kit (50)	50 QIAamp MinElute Columns, Proteinase K, Carrier RNA, Buffers, Collection Tubes (2 ml)	56504
EZ1® DNA Investigator Kit (48)	Reagent Cartridges, Disposable Filter-Tips, Disposable Tip-Holders, Sample Tubes (2 ml), Elution Tubes (1.5 ml), Buffer G2, Proteinase K, Carrier RNA	952034
MinElute PCR Purification Kit (50)*	50 MinElute Spin Columns, Buffers, Collection Tubes (2 ml)	28004

For up-to-date licensing information and product-specific disclaimers, see the respective QIAGEN kit handbook or user manual. QIAGEN kit handbooks and user manuals are available at www.qiagen.com or can be requested from QIAGEN Technical Services or your local distributor.

* Larger kit sizes available; please inquire.

Notes

Trademarks: QIAGEN®, QIAamp®, EZ1®, Investigator®, MinElute®, Rotor-Gene® (QIAGEN Group); 3500™, 6-FAM™, ABI PRISM®, AmpFESTR®, Applied Biosystems®, Avant™, GeneAmp®, GeneMapper®, GeneScan®, Genotyper®, Hi-Di™, Minifiler™ POP-4™, SEfiler Plus™ (Applied Biosystems Corporation or its subsidiaries). Eppendorf®, Mastercycler® (Eppendorf AG); GenBank® (US Department of Health and Human Services); Powerplex® (Promega Corporation). Registered names, trademarks, etc. used in this document, even when not specifically stated as such, are not to be considered unprotected by law.

Notice to purchaser: Limited License

Use of this product is covered by one or more of the following US patents and corresponding patent claims outside the US: 5,789,224, 5,618,711, 6,127,155 and claims outside the US corresponding to expired US Patent No. 5,079,352. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser's own internal research. No right under any other patent claim, no right to perform any patented method, and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel. This product is for research use only. Diagnostic uses under Roche patents require a separate license from Roche. Further information on purchasing licenses may be obtained by contacting the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.

Limited License Agreement for the Investigator DIPplex Kit

Use of this product signifies the agreement of any purchaser or user of the product to the following terms:

1. The product may be used solely in accordance with the protocols provided with the product and this handbook and for use with components contained in the kit only. QIAGEN grants no license under any of its intellectual property to use or incorporate the enclosed components of this kit with any components not included within this kit except as described in the protocols provided with the product, this handbook, and additional protocols available at www.qiagen.com. Some of these additional protocols have been provided by QIAGEN users for QIAGEN users. These protocols have not been thoroughly tested or optimized by QIAGEN. QIAGEN neither guarantees them nor warrants that they do not infringe the rights of third-parties.
2. Other than expressly stated licenses, QIAGEN makes no warranty that this kit and/or its use(s) do not infringe the rights of third-parties.
3. This kit and its components are licensed for one-time use and may not be reused, refurbished, or resold.
4. QIAGEN specifically disclaims any other licenses, expressed or implied other than those expressly stated.
5. The purchaser and user of the kit agree not to take or permit anyone else to take any steps that could lead to or facilitate any acts prohibited above. QIAGEN may enforce the prohibitions of this Limited License Agreement in any Court, and shall recover all its investigative and Court costs, including attorney fees, in any action to enforce this Limited License Agreement or any of its intellectual property rights relating to the kit and/or its components.

For updated license terms, see www.qiagen.com.

© 2011–2014 QIAGEN, all rights reserved.

www.qiagen.com

Australia ■ techservice-au@qiagen.com

Austria ■ techservice-at@qiagen.com

Belgium ■ techservice-bnl@qiagen.com

Brazil ■ suportetecnico.brasil@qiagen.com

Canada ■ techservice-ca@qiagen.com

China ■ techservice-cn@qiagen.com

Denmark ■ techservice-nordic@qiagen.com

Finland ■ techservice-nordic@qiagen.com

France ■ techservice-fr@qiagen.com

Germany ■ techservice-de@qiagen.com

Hong Kong ■ techservice-hk@qiagen.com

India ■ techservice-india@qiagen.com

Ireland ■ techservice-uk@qiagen.com

Italy ■ techservice-it@qiagen.com

Japan ■ techservice-jp@qiagen.com

Korea (South) ■ techservice-kr@qiagen.com

Luxembourg ■ techservice-bnl@qiagen.com

Mexico ■ techservice-mx@qiagen.com

The Netherlands ■ techservice-bnl@qiagen.com

Norway ■ techservice-nordic@qiagen.com

Singapore ■ techservice-sg@qiagen.com

Sweden ■ techservice-nordic@qiagen.com

Switzerland ■ techservice-ch@qiagen.com

UK ■ techservice-uk@qiagen.com

USA ■ techservice-us@qiagen.com

