

Spectralink 87-Series Wireless duet Telephone

User Guide

PIVOT by Spectralink



Copyright Notice

© 2014-2015 Spectralink Corporation All rights reserved. Spectralink™, the Spectralink logo and the names and marks associated with Spectralink's products are trademarks and/or service marks of Spectralink Corporation and are common law marks in the United States and various other countries. All other trademarks are property of their respective owners. No portion hereof may be reproduced or transmitted in any form or by any means, for any purpose other than the recipient's personal use, without the express written permission of Spectralink.

All rights reserved under the International and pan-American Copyright Conventions. No part of this manual, or the software described herein, may be reproduced or transmitted in any form or by any means, or translated into another language or format, in whole or in part, without the express written permission of Spectralink Corporation.

Do not remove (or allow any third party to remove) any product identification, copyright or other notices.

Notice

Spectralink Corporation has prepared this document for use by Spectralink personnel and customers. The drawings and specifications contained herein are the property of Spectralink and shall be neither reproduced in whole or in part without the prior written approval of Spectralink, nor be implied to grant any license to make, use, or sell equipment manufactured in accordance herewith.

Spectralink reserves the right to make changes in specifications and other information contained in this document without prior notice, and the reader should in all cases consult Spectralink to determine whether any such changes have been made.

NO REPRESENTATION OR OTHER AFFIRMATION OF FACT CONTAINED IN THIS DOCUMENT INCLUDING BUT NOT LIMITED TO STATEMENTS REGARDING CAPACITY, RESPONSE-TIME PERFORMANCE, SUITABILITY FOR USE, OR PERFORMANCE OF PRODUCTS DESCRIBED HEREIN SHALL BE DEEMED TO BE A WARRANTY BY SPECTRALINK FOR ANY PURPOSE, OR GIVE RISE TO ANY LIABILITY OF SPECTRALINK WHATSOEVER.

Warranty

The *Product Warranty and Software License and Warranty* and other support documents are available at <http://support.spectralink.com>.

Contact Information

US Location

800-775-5330

Spectralink Corporation
2560 55th Street
Boulder, CO 80301

info@spectralink.com

European Location

+45 7560 2850

Spectralink Europe ApS
Langmarksvej 34
8700 Horsens, Denmark

infodk@spectralink.com

Contents

Chapter 1: Understanding Your Handset	7
<i>Carrying Accessories</i>	7
<i>Handset range</i>	7
PIVOT Hardware.....	8
<i>Front</i>	8
<i>Back 8742 8741/8753</i>	9
Chapter 2: Before Using Your Handset	10
Charging the Battery Pack for the First Use	10
Powering Up and Powering Down your Handset.....	10
Turning on your Handset	11
Turning your Handset Off.....	13
Chapter 3: Getting Started.....	14
Getting Familiar with the Home Screen.....	14
<i>Status bar icons</i>	15
<i>Notifications</i>	15
<i>The Spectralink Dialer</i>	15
<i>Email</i>	16
<i>Apps and Widgets</i>	17
<i>Calendar</i>	18
<i>Browser</i>	18
Setting Screen Security.....	19
Hearing Aid Compatibility	21
Screen Behavior	21
<i>Proximity sensor</i>	21
<i>Controls</i>	22
<i>Auto-rotate</i>	22
Registration Information	22
Chapter 4: Using Telephony Applications	24
Making an Emergency Call.....	25
Ways to Find a Number to Call	26
Call Management	31
<i>Incoming calls</i>	31
<i>Active call</i>	32
<i>Multiple active calls</i>	33
<i>Transfer a call (Blind)</i>	35
<i>Transfer a call (Consultative)</i>	36

<i>Conference calls</i>	36
<i>Call forwarding</i>	38
<i>Speed Dialing</i>	40
<i>Retrieving voicemail</i>	41
Using a Headset or the Speakerphone	42
Using your PIVOT while you are in a call	42
Chapter 5: Using the Barcode Reader (8753 only)	43
General Operating Instructions	44
Sending Barcode Data to PC Applications using QBC	45
<i>Barcode scanning in single-endpoint mode</i>	46
<i>Barcode scanning in multi-endpoint mode</i>	46
Chapter 6: Applications	48
<i>Spectralink widget box</i>	48
<i>All Applications Launcher</i>	48
Chapter 7: SAFE and Panic Button	50
How the SAFE Application Works	50
<i>Icons and notifications</i>	51
<i>Suspending the motion sensor</i>	53
<i>The Warning state</i>	53
<i>The Alarm state</i>	54
Using the Panic Button	55
<i>Panic Button icon</i>	55
<i>Activating the Panic Button</i>	56
Chapter 8: Push-to-Talk (PTT)	58
What is PTT?	58
Using the Push-to-Talk Feature	59
<i>The PTT Application</i>	59
<i>What is a PTT broadcast?</i>	61
<i>Transmitting a PTT broadcast</i>	61
<i>Receiving PTT Broadcasts</i>	64
<i>Receiving PTT Broadcasts while in Call</i>	66
<i>Receiving Priority PTT Broadcasts</i>	67
<i>Receiving Emergency PTT Broadcasts</i>	68
<i>Changing the Default PTT Channel</i>	68
Chapter 9: Chargers and Battery Packs	70
<i>How usage affects Battery Pack life</i>	70
Charging the Battery Pack	71

Charger Options	72
<i>Using the USB Charger.....</i>	<i>72</i>
<i>Using the Quad Charger</i>	<i>73</i>
About Battery Packs.....	74
<i>Indications of Low Battery</i>	<i>74</i>
<i>Battery Pack Removal and Replacement</i>	<i>75</i>
Battery Safety Notices.....	76
<i>Notes on Battery Packs.....</i>	<i>76</i>
Chapter 10: Sync Calendar and Email	78
Chapter 11: Settings	83
Settings Options	83
<i>Sound</i>	<i>83</i>
<i>Display</i>	<i>84</i>
<i>Battery</i>	<i>84</i>
<i>Language & input.....</i>	<i>85</i>
<i>Add an account</i>	<i>86</i>
<i>Accessibility</i>	<i>86</i>
<i>About phone.....</i>	<i>87</i>
Updating Handset Software	88
<i>Verifying the update</i>	<i>88</i>
Chapter 12: Other Things You Should Know	89
<i>No camera?</i>	<i>89</i>
<i>Can't access the Web?</i>	<i>89</i>
<i>Can't connect to a computer?.....</i>	<i>89</i>
<i>Can't print?.....</i>	<i>89</i>
Appendix A: Care of the Handset	90
Appendix B: Regulatory and Safety Information	92
Product Compatibility.....	92
FCC	93
Industry Canada Notice	94
CE Declaration of Conformity	95
Hearing Aid Compatibility	98
Chargers and Battery Packs	98
Wireless Telephones	99
<i>Operational Warnings</i>	<i>100</i>
Electromagnetic Interference and Compatibility	100
<i>Facilities.....</i>	<i>100</i>
<i>Pacemakers.....</i>	<i>100</i>

<i>Other Medical Devices</i>	101
<i>Use While Driving</i>	101
<i>FCC Statement</i>	101
Specific Absorption Rate (SAR) Information: Spectralink 87-Series Wireless	
Handsets	101
<i>EU SAR peak values (in W/kg):</i>	102
<i>FCC SAR peak values (in W/kg):</i>	102
<i>Handset operation normal position</i>	103
Appendix C: Products Mentioned in this Document	104

Chapter 1: Understanding Your Handset

Your PIVOT handset combines the familiarity of the Android™ platform and Google™ applications with the advanced functionality and durability needed for a workplace environment. PIVOT is not a cell phone. It is a Wi-Fi device that utilizes the wireless infrastructure installed in your workplace for connection and communication. We call it a WorkSmart device.

The PIVOT family of products consists of two form factors and three models. All models feature a large, backlit color display, a touchscreen with full keyboard, an enterprise-grade speakerphone and superior resistance to dust and liquids, including cleaning agents. The PIVOT:S (8742) has a streamlined design with a slim profile for a familiar smartphone feel. PIVOT 8741 has a patented kickstand spine design. PIVOT 8753 includes an integrated, high-performance 1D/2D barcode scanner. All models have removable, rechargeable batteries.

Like other devices powered by Android, PIVOT uses applications to deploy functionality. Each facility chooses what applications it makes available to its PIVOT users. Contact your system administrator for information about what capabilities are configured for your facility.

This user guide shows you how to operate your Spectralink 87-Series Wireless Handset. Additionally, this chapter orients you to the physical features of the handset.

Spectralink applications are translated into French, German, Spanish and Italian but we have retained the long list of languages offered by the Android platform. If you select a different language, labels associated with specialized applications may appear in English.

Instructions for using familiar Android features are readily available online. Specialized Spectralink applications follow these same usage guidelines. Use this document as a starting place to explore PIVOT and learn a new way of working smart.

Carrying Accessories

PIVOT:S 8742

- Use a holster to carry and protect your PIVOT:S model.

PIVOT 8741 and 8753

- Two clip styles are available: Both the clothespin-style and button-style easily attach to a belt or other firm edge.
- Cross-body lanyard hugs the body for secure carrying and quick access.
- Silicone bumper protects your PIVOT without interfering with use.

Handset range

Your handset will work within your facility only when it is in contact with the WLAN and Local Area Network (LAN). It will not work beyond the RF coverage area of the wireless LAN. Contact your system administrator if you have questions about your facility's RF coverage area.




Note: Configurable features

Many of the features detailed in this Guide are configurable by your system administrator and may not follow the exact steps outlined here. Contact your system administrator for more information if you have trouble understanding how to use certain features.

PIVOT Hardware


Front

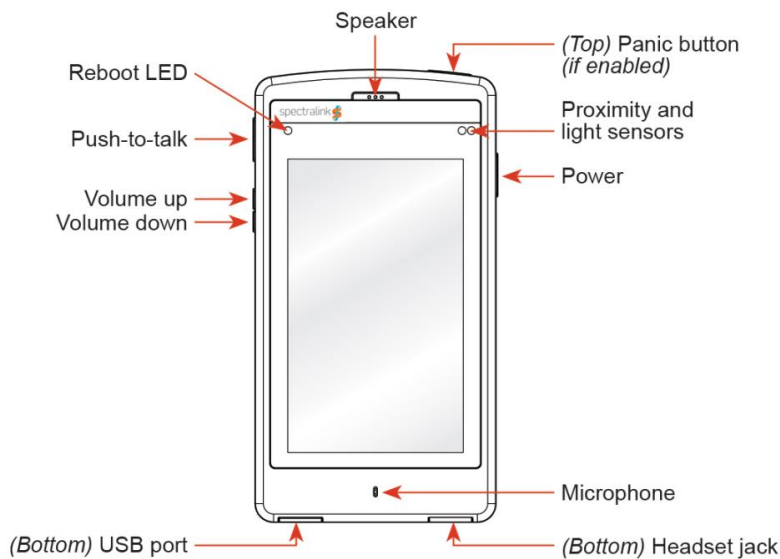


Note: The Push-to-talk and Panic buttons

The Push-to-talk and Panic buttons each require a separate application that may or may not be enabled in your facility. Consult with your system administrator for complete information.

Back 8742

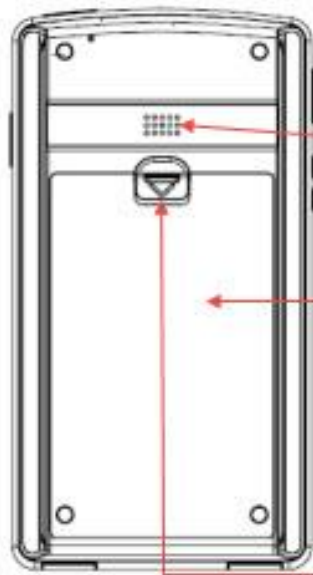

8741/8753

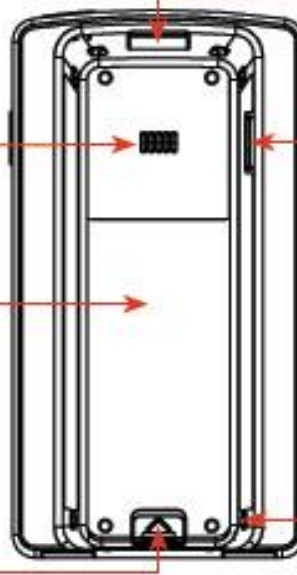

Rear speaker

Scan window (8753 only)

Scan button (8753 only)

Battery Pack

Battery Pack latch

(All 4 corners)
Lanyard through-holes

Chapter 2: Before Using Your Handset

When you first receive your handset, you need to fully charge the Battery Pack in order to maximize the Battery Pack's storage capacity and lifespan. You may also need to familiarize yourself with turning the handset on and off and learn about the various audio options available to you. Additionally, your facility may have security standards that require you to enter your unique credentials so that the handset will correctly register with your system. This chapter covers all of these subjects.



Note: Indoor use only

This device is intended for indoor use only.

Charging the Battery Pack for the First Use

Take a moment to learn how to charge your handset properly. It is critical to fully charge the Battery Pack before the first use of the handset in order to maximize the Battery Pack's storage capacity and lifespan. The Battery Pack is fully charged when the Battery life indicator becomes solid.

Battery is plugged in



Battery is fully charged



Caution: Always charge your handset when it is ON.

Your PIVOT handset must be turned on in order for it to accept a charge. Do not turn your handset off when you plug it in for charging with the USB cable.

Powering Up and Powering Down your Handset

If the Battery Pack has adequate charge, your PIVOT will turn on when the power button is pressed.

To power down the handset, press and hold the power button until the power off popup appears.



Warning: Loss of data

Do not remove the Battery Pack to power off the handset, as you may lose personal data. If you need to replace the Battery Pack, power off the handset and then remove the Battery Pack.

Turning on your Handset

Generally, your system administrator will deploy your handset with multiple other handsets. In this case, your system administrator will configure the handset parameters beforehand so that the handset will register with the system and be ready to use when you turn it on. If your handset does not work when you turn it on, contact your system administrator.

- 1 Press the power button to turn on the handset



Note: LED behavior

When you turn on the handset, the LED indicates normal and abnormal operating conditions. Normally, the LED will come on, the handset will vibrate briefly, the display will turn on, and the LED will turn off. Other circumstances may cause the handset to deviate from this behavior.

LED off and then flashing 5 times, handset does not turn on: The Battery Pack must be charged.

LED off and then flashing 10 times: Check to see if the Battery Pack is properly seated. Replace if flashing behavior continues. Give faulty Battery Pack to your system administrator.


LED off, then on for two seconds then handset turns off: Indicates an internal error. Try to power on the handset again. If problem persists, contact your system administrator.

LED on, off briefly, then on again several times: The battery is being checked.


- 2 The Spectralink logo displays followed by the PIVOT logo.
- 3 The lock screen displays showing the lock icon that the system administrator has configured. The default lock method is swipe.

Place your finger on the lock icon and swipe towards the unlock icon shown to the right of the lock icon.

Lock screen showing the slide lock

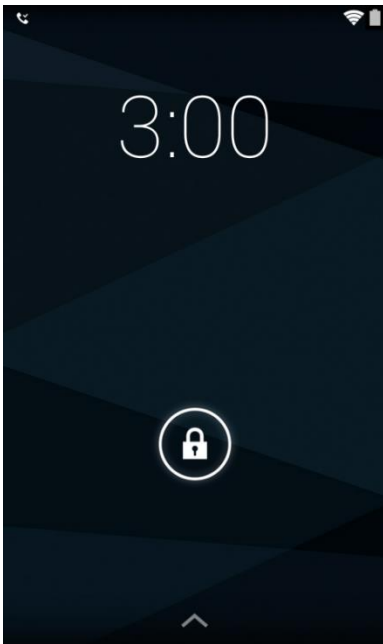

Lock screen showing the PIN lock

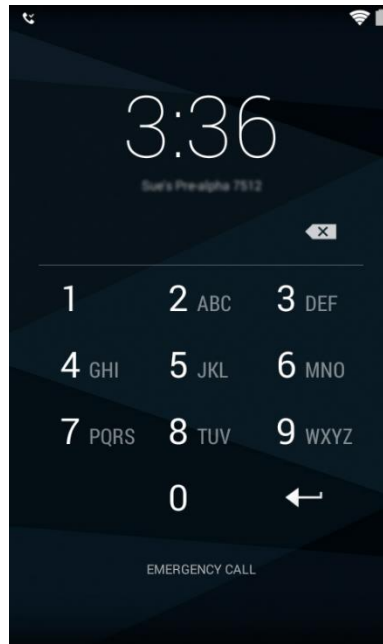

Note: Emergency call safety feature


If your handset is configured with a PIN, Pattern or Password lock, the Emergency call button appears, enabling you to make a call before unlocking the handset. Tap the button to open a list of emergency contacts if they have been configured by the system administrator.



Note: User information

If configured by the administrator, user information will display below the date.

- 4 If credentials have been configured, you must enter your PIN or Password before you can use the handset. Simply tap the keypad/keyboard to enter the characters.
- 5 The home screen displays.

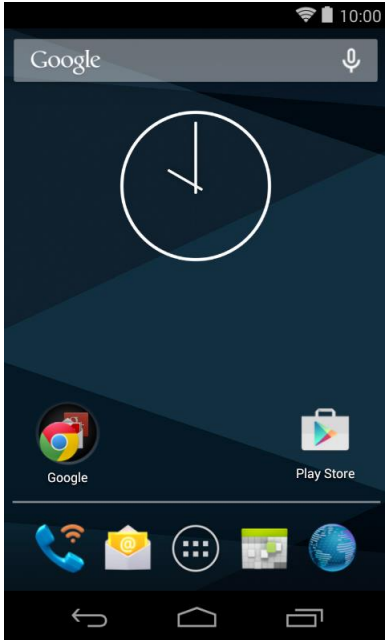


Your home screen may differ depending on the customizations done for your facility.

- 6 Pressing the power button briefly will turn the screen off and lock the handset. Press the power button briefly again to turn the screen back on.


Turning your Handset Off

To turn your handset off

- 1 Press and hold the power button. The Power off popup appears:

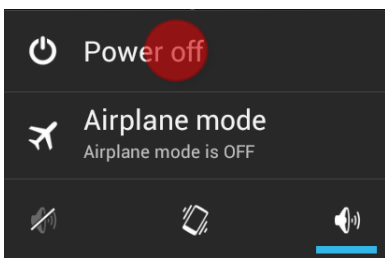

- 2 Tap Power off. The confirmation popup appears:

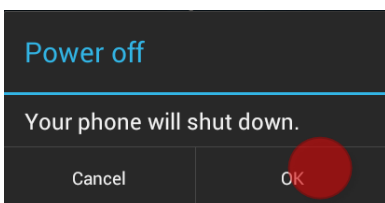

- 3 Confirm by tapping OK

Chapter 3: Getting Started

This chapter shows you how to use the Home screen, the Status Bar icons and generally how to navigate through the PIVOT screens.

Like other devices powered by Android, PIVOT is app driven, not menu driven. Generally you will tap icons to open applications that provide features, such as dialing a call, viewing email, browsing the Web, etc.




Note: Settings icon

Like other devices powered by Android, you can customize many settings, such as phone ring tones and volumes, through a Settings application which is activated at the tap of its icon. You can find the Settings app on one of your home screens and in the Launcher screen.

Getting Familiar with the Home Screen

The Home screen is the first display you will see. At the top, the Status Bar displays icons that give you information about the status of your handset.

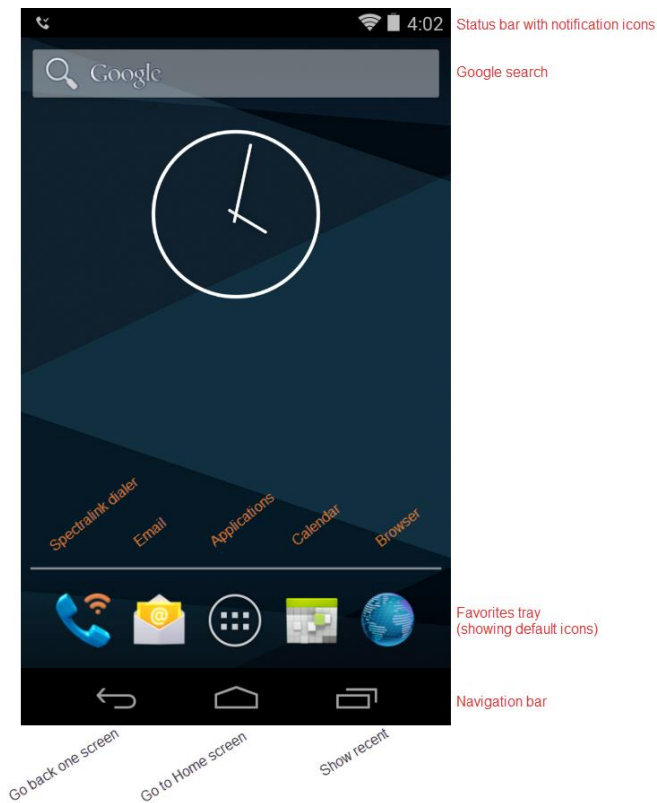

Status bar icons

The status bar is dynamically populated to provide current information about your handset. Minimally you will see the Wi-Fi icon, the battery icon and the current time.



Your Pivot is a Wi-Fi device that uses a wireless connection to the network in your facility. The Wi-Fi icon indicates signal strength of the connection. Low signal strength will prevent proper handset operation. Contact your system administrator if you do not have sufficient signal strength as indicated by grayed out or no bars on this icon.



The battery icon represents the amount of charge remaining in the Battery Pack.




The digital clock draws its time from an online resource or from a local time server. See your system administrator if the displayed time is incorrect.

Notifications

Notifications appear on the left side of the status bar. Missed call and voicemail notifications are an example. The handset icon indicates that the handset is registered and you can make calls.




You can swipe down to open the notification drawer and get additional information.

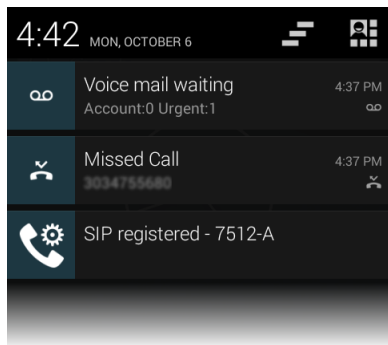

When the handset is idle and the screen is black, a blinking LED lets you know that you have a notification.

The Spectralink Dialer

The Dialer icon opens the Spectralink dialer application. The keypad displays by default.

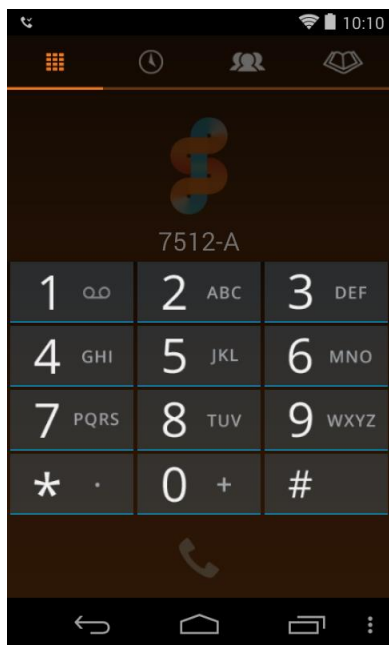

Four icons across the top of the screen open different tabs:



The keypad tab opens the keypad. It is the default tab.



The recent tab opens a list of recent calls.



The contacts tab opens a list of contacts as populated by the People app.




The corporate directory tab opens the directory of corporate contacts if activated by you system administrator.

Go to [Using Telephony Applications](#) for complete information.

Email



The email icon allows you to setup an email account for sending and receiving email messages. Contact your system administrator for information about setting up your corporate email account. See [Sync Calendar and Email](#).

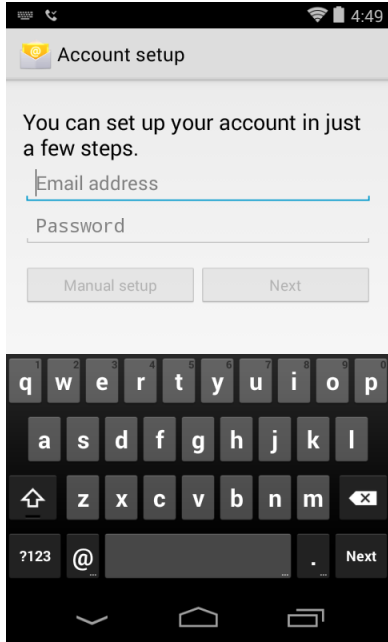


Apps and Widgets




Like other devices powered by Android, you can customize your home screens with apps and widgets by dragging them to a Home screen from the Launchers.

Apps Launcher

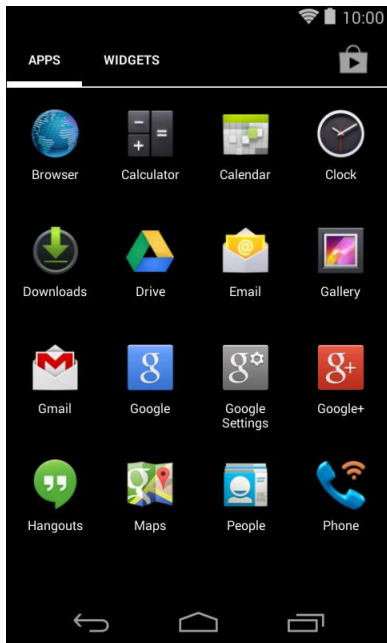


Widgets Launcher

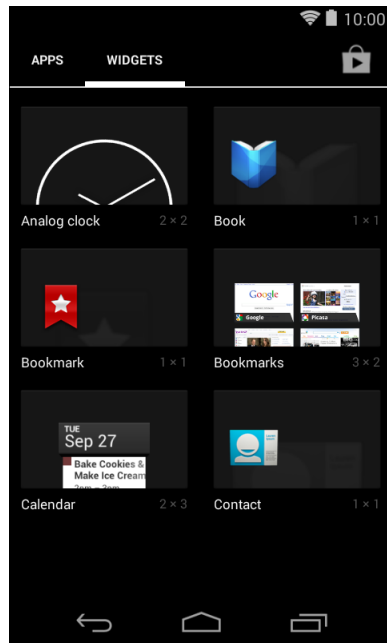


These can be added to and removed from any Home screen as you would with any device powered by Android: To add an app, press and hold the icon and then drag it when the Home screens appear. To remove an app from a Home screen, press and hold the icon and then drag it to the **X Remove** area when it appears.

Calendar

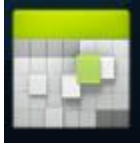


The Calendar can be synced with your Outlook program. Contact your system administrator for information on syncing Outlook. Set up is described in [Sync Calendar and Email](#).

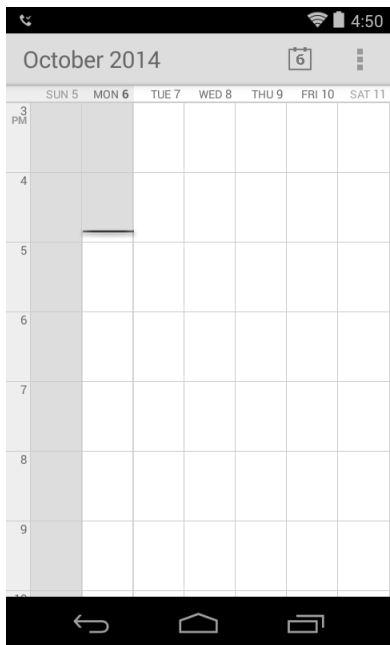

Browser



The Browser icon opens the search page on Google (or other home page as set by your system administrator).

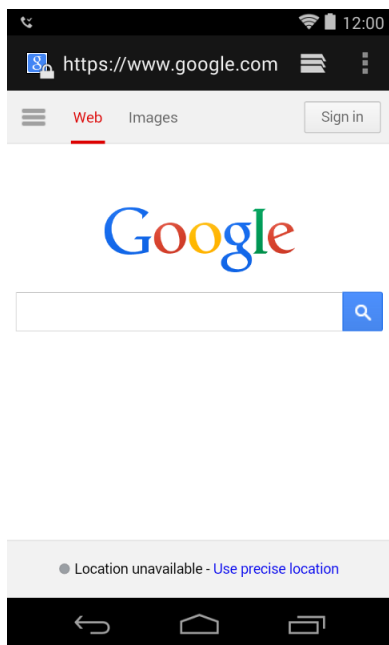

Setting Screen Security


To avoid making inadvertent calls and to prevent unauthorized calls PIVOT offers several lock methods. When activated you must unlock the keypad before you can make calls. Your system administrator may have set a lock method already or may give you the option of setting a secure method. Be aware that if you go from an insecure method—Slide or None—to a secure method—Pattern, PIN or Password—you will see a warning that returning to an insecure method requires help from your system administrator. Check with your administrator if you need more information about screen lock requirements in your facility.

Some installations require that you change your Pattern, PIN or Password periodically. Use these directions to make this type of change.


Check with your system administrator if the Screen security feature is not available and for information on the timeout value.

If your system administrator allows you to change the lock setting, navigate to **Settings> Screen security** and tap **Screen lock**.


Settings menu

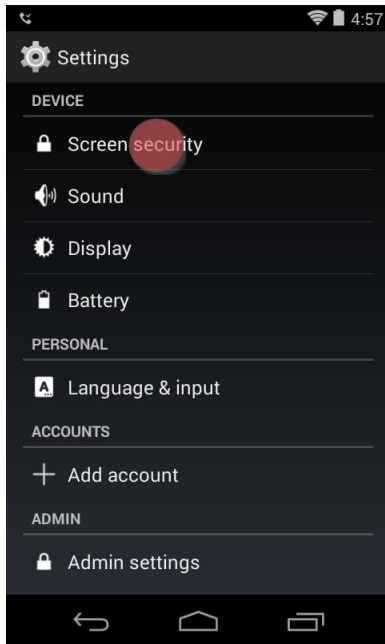

Screen security option

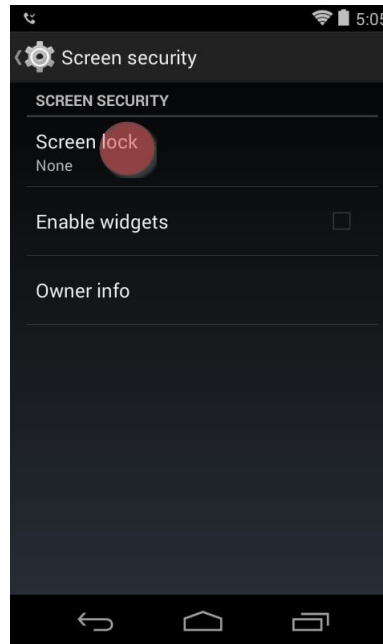

Choose screen lock

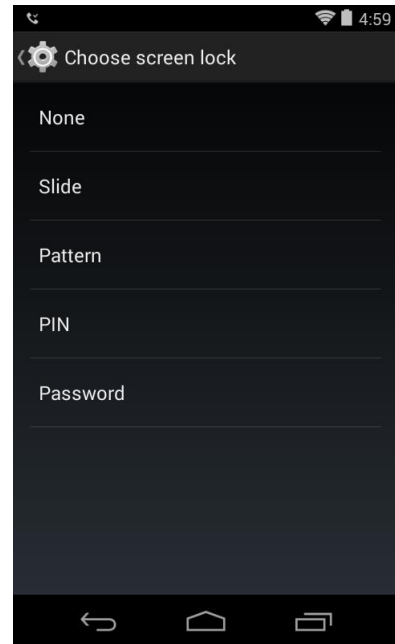


Note: When lock options are disabled

Your system administrator may configure your handset with a Microsoft® Exchange Active Sync email account which requires you to use a PIN or Password to enable handset operation. If so, the None, Slide and Pattern options will not be available.

When you change the lock, additional options will be shown:

PIN screen lock options

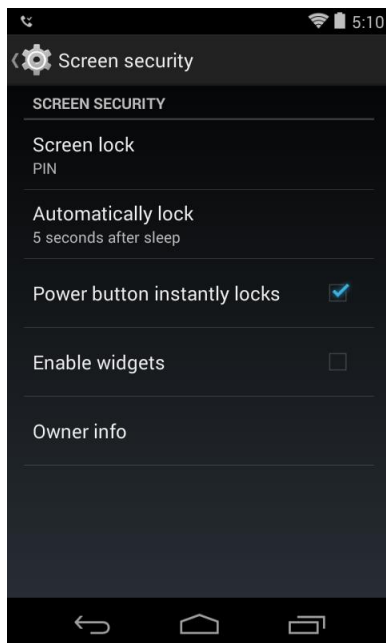

What is that clicking sound?

If the system administrator has activated encryption on your handset, you may hear a clicking sound until you enter your unlock credentials the first time after powering up your handset.

Hearing Aid Compatibility

Your handset has a hearing aid compatibility mode that adjusts the frequency response for users who have hearing aids equipped with telecoil inductive pickups. Acoustic pickup hearing aids will not respond to this adjustment.

To view the current status of the Hearing Aid Compatibility feature and to enable or disable this feature, navigate to **Settings> Sound> Hearing Aid Compatibility**. Tap to activate. A checkmark indicates that the option is activated.

Screen Behavior

Proximity sensor

The proximity sensor is located on the top right corner of the handset. It senses when the handset is close to your ear and turns off the screen so that you do not inadvertently tap a button with your cheek. When you move the handset away from your head, the screen turns back on and will accept your finger taps.

Controls


The standard Android platform used by PIVOT offers finger motions that are used to manage the screen such as tap, swipe, long press, press and hold, drag, etc., depending on the situation.


Auto-rotate


You can change the screen to landscape display during some applications. Toggle this capability on and off through the Settings menu. See [Display](#).


Registration Information

Your PIVOT is dependent upon both the wireless and wired network in your facility. When it is registered to a call server and the network is working properly, you can make and receive calls and perform various other telephony functions as programmed in your facility.

The registration icon  displays in the status bar whenever registration is successful. A notification will display in the notification drawer:

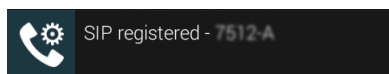


If registration fails, you will have the failed registration icon in the status bar  and a notification in the notification drawer.

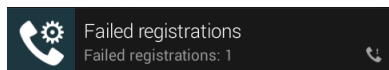

Check registration information by going to the Launcher screen and tapping the SIP Status icon.



Successful registration

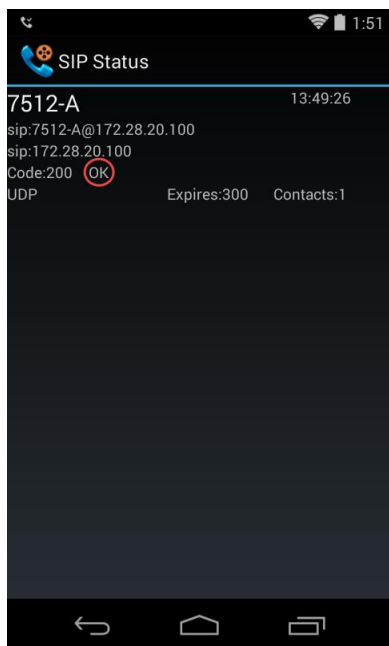

The SIP Status screen gives you call server registration information. Your system administrator may ask you for this information when you call for support.

Chapter 4: Using Telephony Applications



Warning: Power failure may prevent call functionality


Your handset may not function in the event of a power failure as it relies on the greater computer network to transmit calls. The Wi-Fi icon will be completely grayed out if the wireless system is down. You will have a registration failure icon if the handset cannot make or receive calls, but there may be a time delay between the system failure and the appearance of the registration failure icon.

Check with your system administrator to see if the wireless telephony system is protected from power outages.

The Spectralink Dialer provides telephone call functionality. Tap the Dialer icon on the home screen to open the Dialer.



The action bar across the top displays the action buttons you can use for direct dialing or to select people to call.

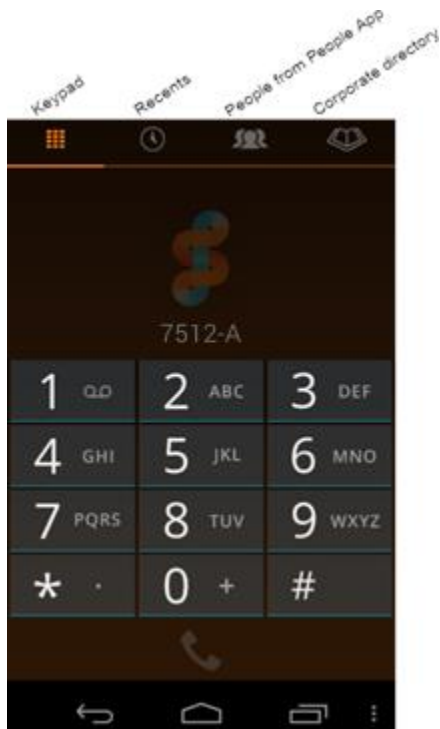


Making an Emergency Call



Settings: Emergency calling options while the keypad is locked

If a lock option requires a pattern, PIN or password you can make an emergency call if emergency numbers have been configured by the system administrator.

Tap Emergency call below the input pad and tap the number from the list of numbers configured for your facility.

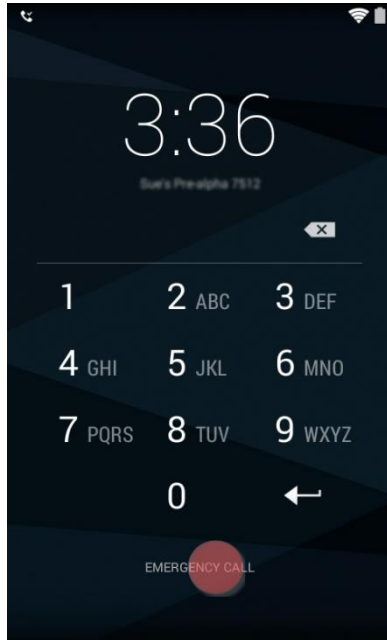

Contact your system administrator for additional information about numbers you should call in an emergency situation and how to access them with your PIVOT.

Ways to Find a Number to Call



Note

There are many different ways to populate the dialer, several of which are detailed below. Once populated, in order to dial the number you must tap the Start button.



Caution: Dialing anomalies

Your facility may have certain dialing requirements that are not reflected in numbers that are captured from call history, the directory or even numbers entered in the People app.

Such required numbers must be entered manually after the Dialer is populated. Just tap the field and enter the numbers you need (or delete the numbers you don't need).

PIVOT offers the "+" sign for international dialing if supported in your facility. Enter + in the dialer with a long press of the 0 key.


It is good practice to verify the number you are calling is correct before tapping the Start button.

Dialpad call


When you tap the numbers on the keypad, the dialer will attempt to find a match in the People app and display the resulting information.

In the screens below, there were no matches to the full number entered on the left. On the right side screen, a match from the People app occurred when the number 7 was entered.

No match

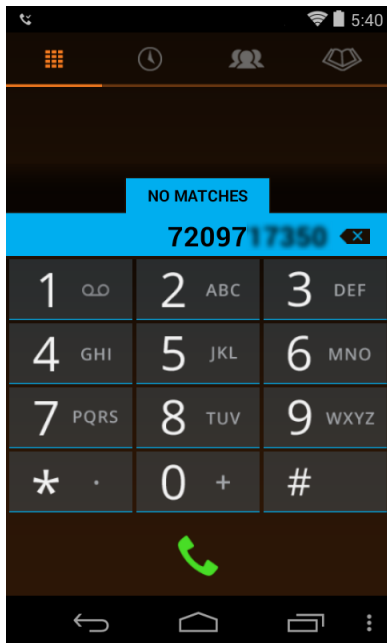


Match from People app

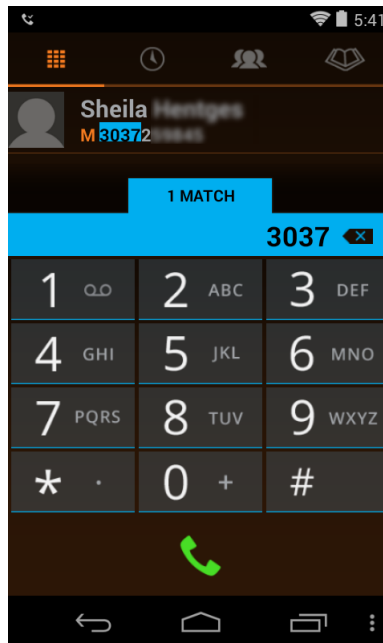

Recent call


Tapping the Recent icon opens a list of recent calls. Tap the one you want to dial. Remember, you may need to add or delete numbers according to the dialing standards in your facility.

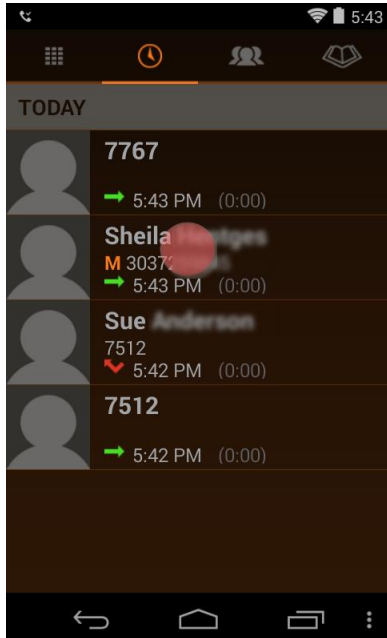



Timesaver


Use a long press on a recent call to add the number to the People app as a contact. When you long press a number and tap Add to Contacts, People opens and you can tap Create New Contact and enter the name and other contact information.

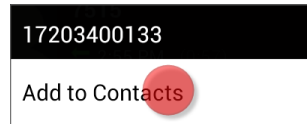


Contacts

When the People app has been populated, you can tap the Contacts icon to open the list of your contacts and tap the contact to populate the dialer. Remember, you may need to add or delete numbers according to the dialing standards in your facility.


Tap the name to populate the dialer

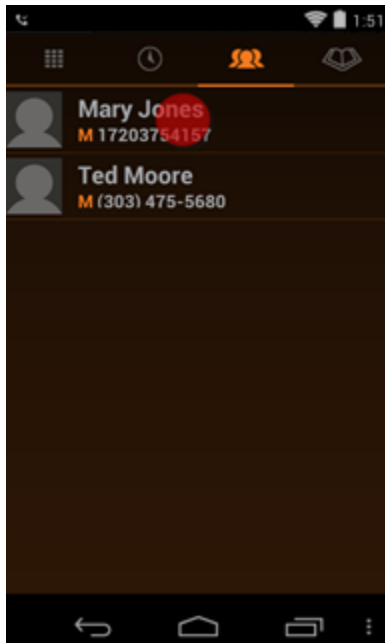

Tap Start to place the call

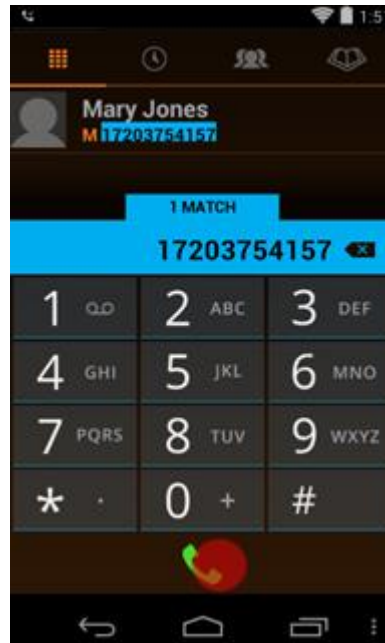

The call is answered...

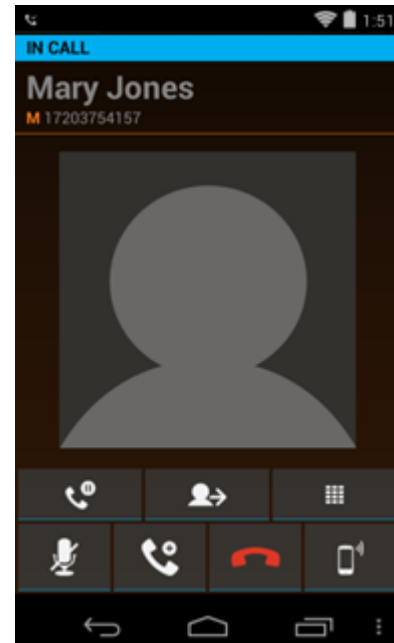


To add contacts to the People app

The People app can be found on the Launcher screen.




- 1 Tap the icon to open the People app and set up your profile.
- 2 Add contacts in three ways: manually entering them, entering them from the Recents list in the dialer, importing them when you sync with Outlook® (see [Sync Calendar and Email](#)).
- 3 You can type in your contacts.


Tap the Create button

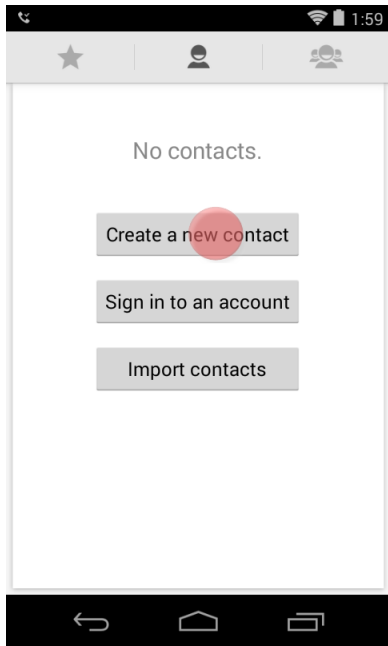

The contact entry

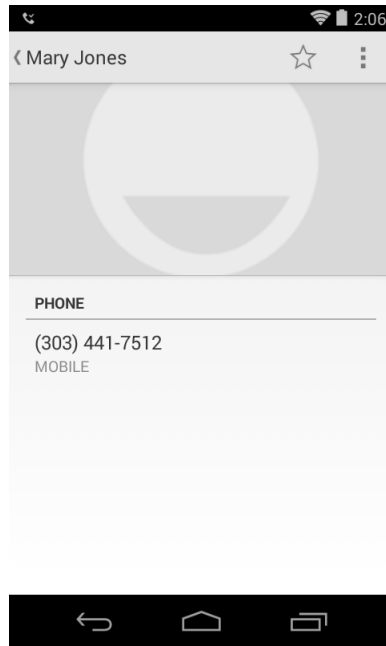

Available from the dialer

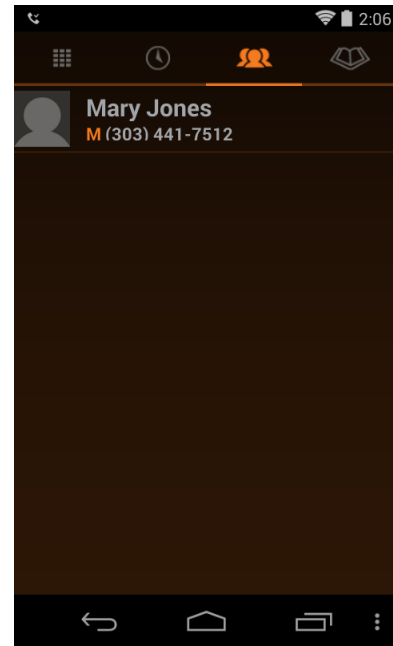

- 4 You can also add contacts from the Recents list by using a long press and tapping Add to Contacts on the popup.




Note: Trouble with People?

If your PIVOT handset has an extremely large number of contacts, you may experience performance issues with the Spectralink phone application. Contact your system administrator for assistance if your contacts list searches slowly or if you get an error message.


Corporate Directory

If installed in your facility, the corporate directory is a searchable database that you can access within the dialer.


Tap to open the Corporate Directory

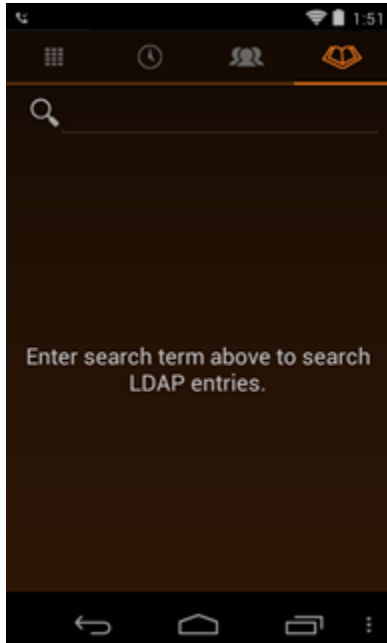


Enter text and tap the search button

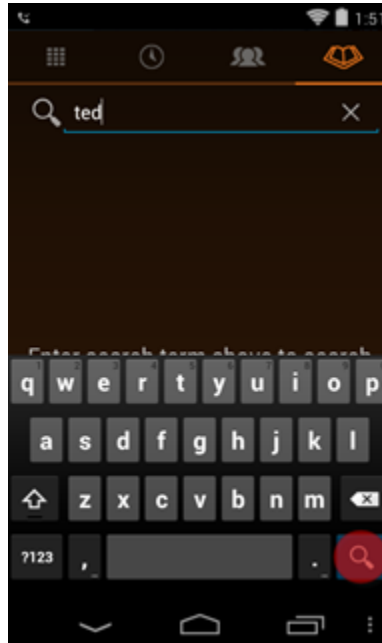

View search results, tap name

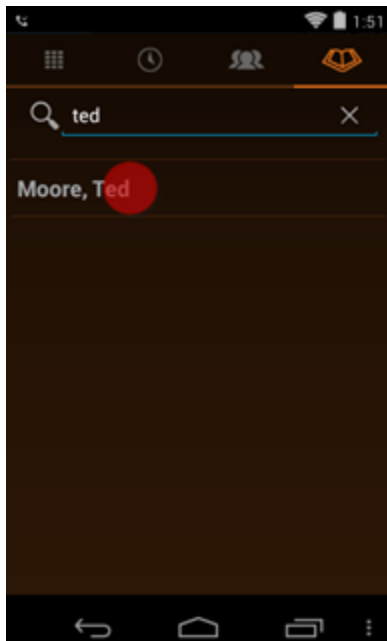

View results, tap number to place call

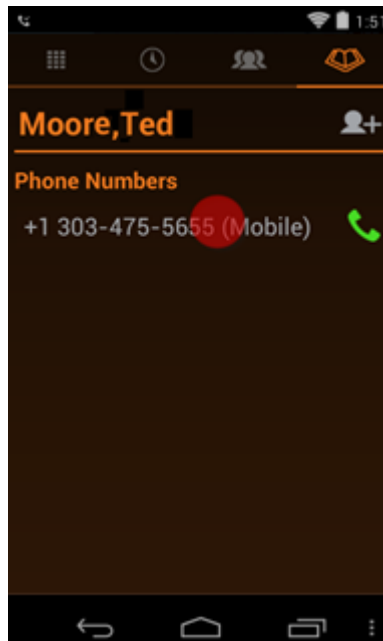


Tapping the phone number opens the dialer.

Call Management

As a device powered by Android, call management on your PIVOT will seem familiar. However, your PIVOT has also been designed to have many features found on your corporate deskphones. As devices powered by Android are slightly different, look through these call management instructions to learn about how PIVOT functions.

Incoming calls

An incoming call screen

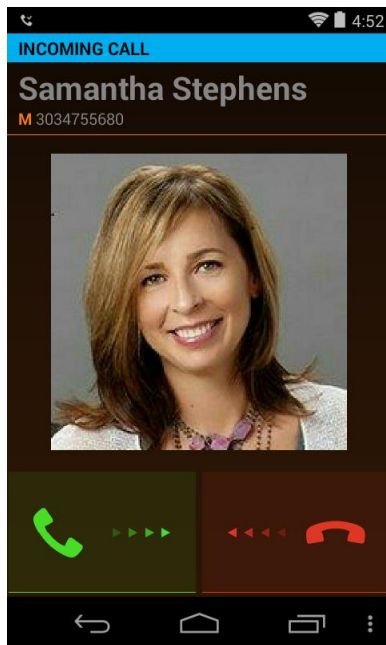

Answer an incoming call

To answer, slide the green button to the right.


Silencing an incoming call

Press a volume button to turn off the incoming ring without ending the call.




Note: Silencing your handset

You can set your handset to silent in the power off popup.

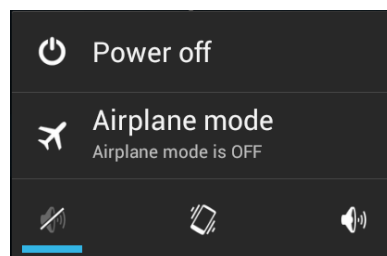


Reject an incoming call

To reject, slide the red button to the left. The call will stop ringing. When you reject a call, the missed call notification will briefly display on the status bar, it will disappear and the missed call icon will remain on the status bar. A voicemail icon will display if a voicemail has been recorded (and if voicemail is activated in your facility).

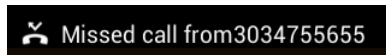

The missed call notification will display in the notification drawer (opened by a down-swipe).

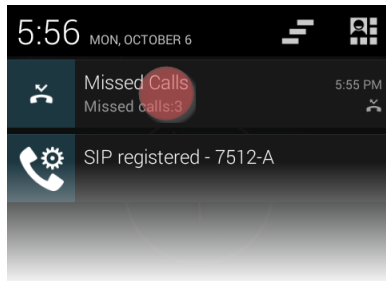

Tap the notification to open the Recents list in the dialer and view the missed call(s). This clears the missed call notification in the status bar.



Note

Your missed call notifications will be erased when the phone is turned off.

Active call


When you answer a call by swiping the green call icon, it is immediately active. Hold the handset to your ear and speak.

An active call screen



During an active call, the action bar displays icons to help you manage a call. Hold, Mute, End and Audio option are fairly self-explanatory. Transfer and Add call (conferencing) are explained below.

The action bar

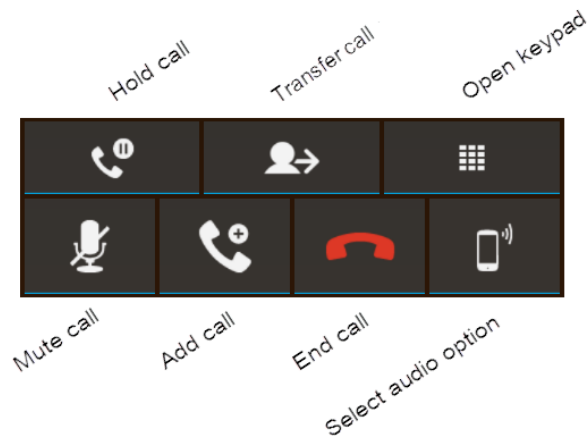


Multiple active calls

If multiple active calls are permitted in your facility, you can have up to four calls at one time. Simply follow steps 1-4. Contact your system administrator if you have questions about how many active calls you can have at once.

Setting up multiple calls

- 1 While in an active call, tap the **Add call** button. This places the first call on hold and opens the Add call dialer.

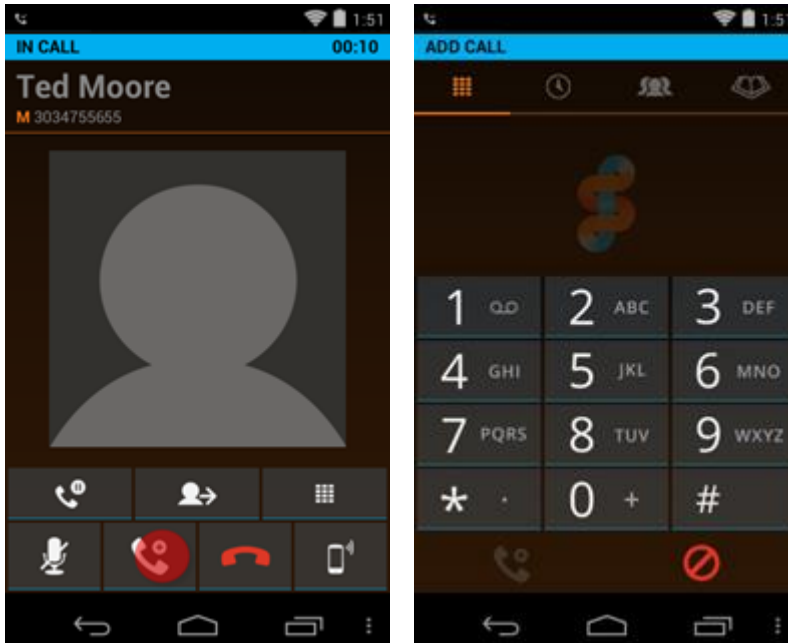


- 2 Use any method to populate the dialer for a call to the target party.



- 3 Tap the green Add call button to dial the number.

- 4 When the target party answers, you will see the active call screen with both calls displayed and the first call on hold. You can alternate between calls by tapping the other call which will place the first call on hold.

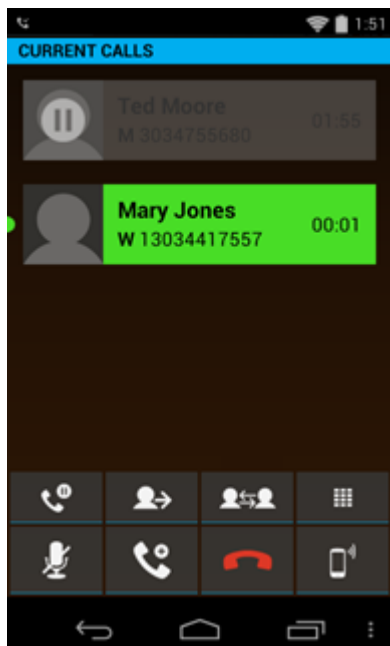


- 5 Follow the above steps to add additional calls. Your PIVOT is capable of being in separate calls with up to four numbers at once.

Transfer a call (Blind)

To transfer a call means to send an active call to another number. When you do this, the call is no longer active on your handset. An unattended or blind transfer is when you do not speak to the party you are transferring the call to before the call is transferred.

Setting up a blind transfer

- 1 While in an active call with the first party, tap the Transfer button. The dialer opens.
- 2 Use any method to populate the dialer.

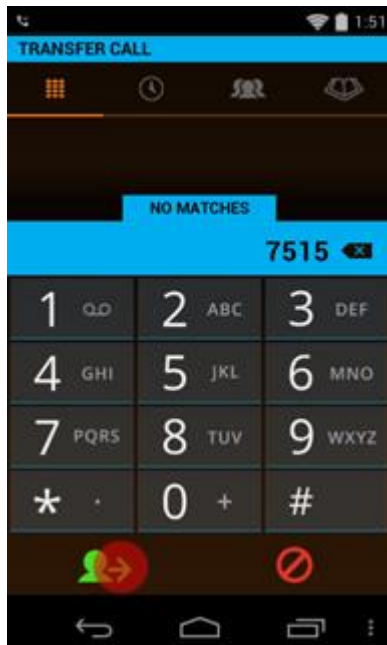

- 3 Tap the Transfer call button to complete the transfer.
- 4 Your screen will return to the Home screen.

Transfer a call (Consultative)

In a consultative transfer, you are able to speak to either of the two other parties before making the transfer.

Setting up a consultative transfer

- 1 When you are in a call with the first party, add a call to the second party. This puts the first call on hold.
- 2 When the second party answers, explain that you are going to transfer the call.
- 3 Tap the Transfer call button and then tap the first call to complete the transfer.
You are returned to the Home screen or, if other calls are active, to the active call screen.

Conference calls

A conference call is when you and two other parties can converse at once. Your PIVOT supports a maximum number of three participants in a conference call—you and two others.

Set up a conference call


- 1 Follow steps 1-4 under Multiple Calls.

- 2 When the target party answers, you will see the active call screen with both calls displayed and the first call on hold. The call to the target party is active and you can explain that you are going to set up a conference call.



- 3 Tap the **Join** softkey. When all three participants are connected to the conference call, the conference screen will look like this:


Conference call screen example

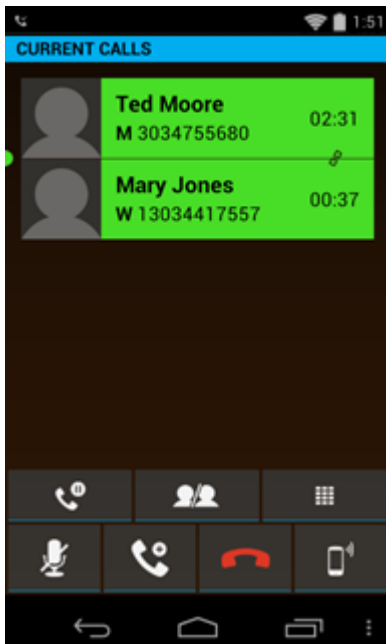


Ending a conference call

- Tap the End call button when both calls are active. This hangs up both called parties.
- Tap the Split button. This splits the conference call into two separate calls, both on hold. You can end the call with one participant and continue the call with the other party. Or you can carry on separate conversations alternately or rejoin them back into a conference call.

Tap the Split button

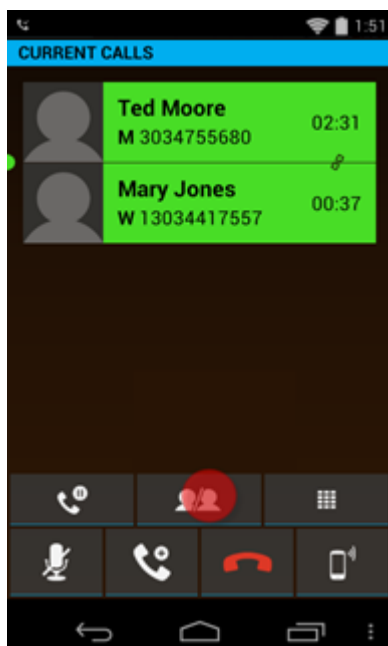

Split conference call showing two single calls on hold

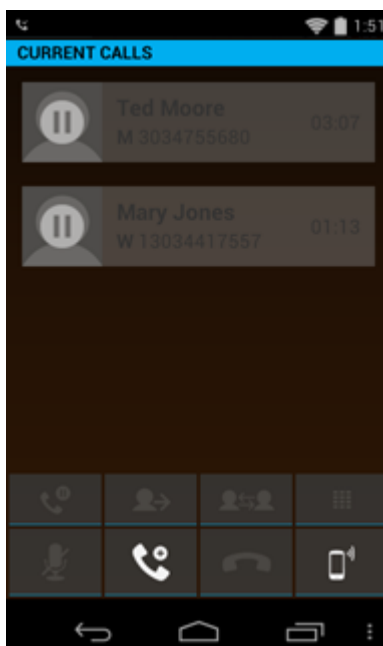


Call forwarding

Call forwarding is when calls to your phone are routed to a different number and your phone does not ring at all.

Set up call forwarding


- 1 Open the dialer and tap the overflow button on the navigation bar. The Call Forwarding Setup button appears in a popup.

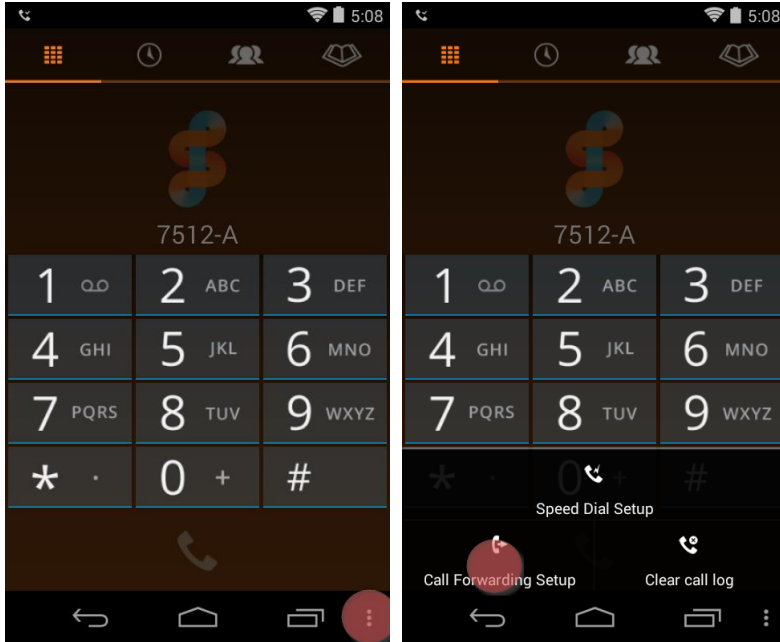

- 2 Tap Call Forwarding Setup.
- 3 Slide the ON/OFF slider to ON and enter a number or tap a contact from the displayed contact list in the People app.

Call Forwarding Setup

Enter number

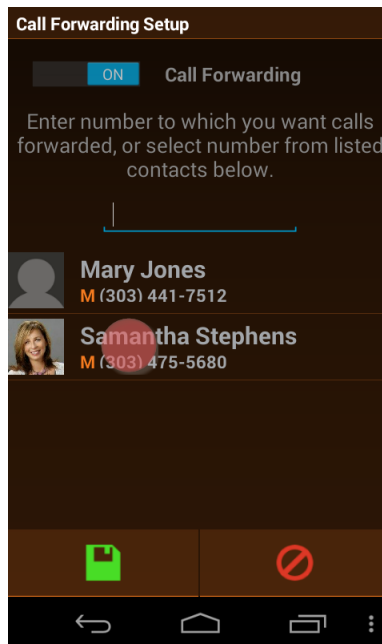

- 4 Tap the Save button to confirm the number. Call forwarding is now activated.

Turn off call forwarding

Tap the overflow button and open Call Forwarding Setup. Slide the ON/OFF slider to OFF and tap Save. Call Forwarding is now disabled.



User Tip: Call forwarding

Automatic call forwarding works only as long as your handset is registered and turned on. If you turn off your handset or move out of range for an extended period of time, calls to your handset are neither received nor forwarded. Depending on your organization's Private Branch Exchange (PBX) programming, the caller may be sent to the Voicemail system, hear unending ring tones, or a message such as "The party you have called is not available".

Speed Dialing

When you have a lot of names and numbers in your People list, the speed dial feature is very useful as you can assign a speed dial number to a contact and quickly dial the number. Nine numbers are provided, including number one which may already be assigned to voicemail.




Power User

You must add numbers you wish to use as Contacts in the People App before attempting to program them as Speed Dial entries.


Set up speed dialing

- 1 Open the dialer and tap the overflow button on the navigation bar.

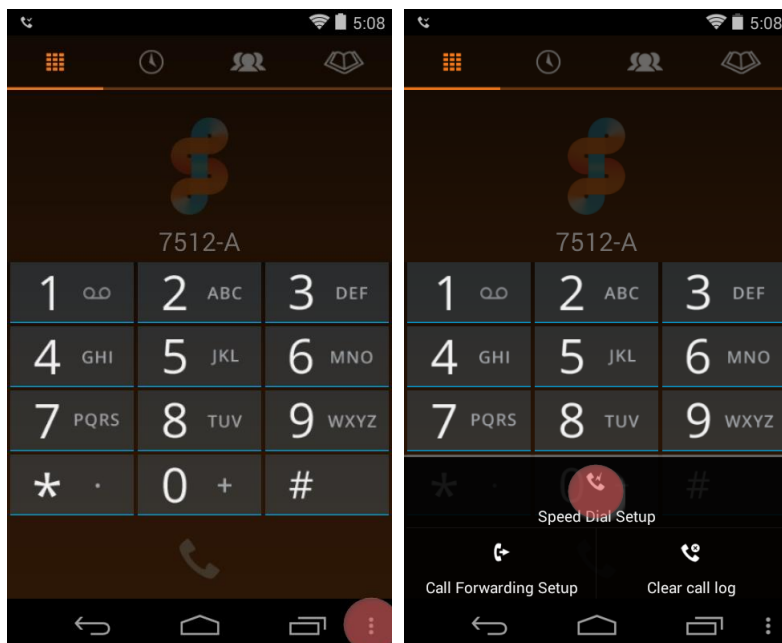


- 2 Tap Speed Dial Setup.
- 3 Tap a speed dial number to assign. The People list opens. Tap the contact number you want to assign to the speed dial number.


Speed Dial Setup

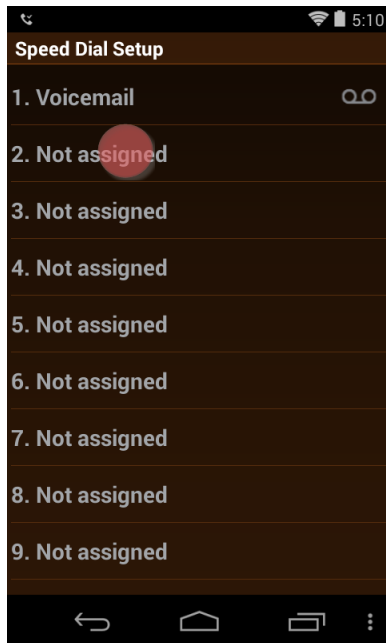

Tap number to assign

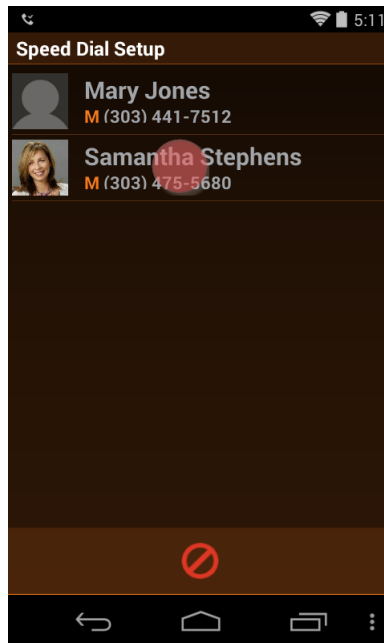

Speed dial assigned

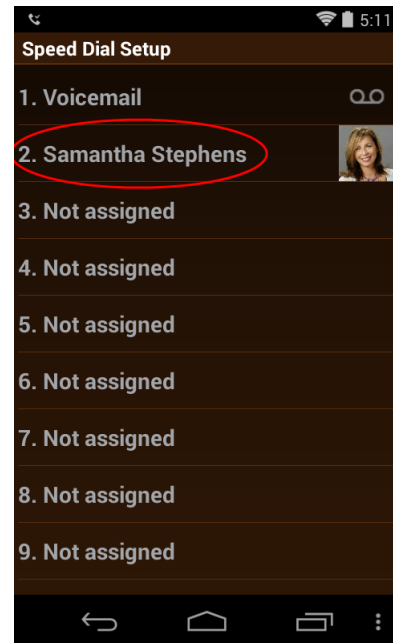

Using speed dial

- 1 Open the dialer keypad.
- 2 Long-press the speed dial number. The call is immediately placed.

Retrieving voicemail


Voicemail may be set up in your facility.

Listening to a new voicemail

A new voicemail will be indicated by the voicemail icon in the status bar.



Swipe down to open the notification drawer. Tap the voicemail waiting notification to call your voicemail retrieval number.

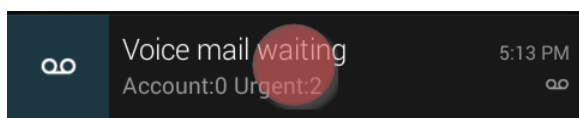



Note

If you hang up during the call to your voicemail retrieval number, the notification will be gone from the status bar, but you will still have a voicemail that you have not yet listened to. In this case, you can dial your voicemail number to retrieve the message.

If your voicemail messages do not show a notification when they should, powering off the phone and then turning it back on will restore notifications.

Call your voicemail retrieval number

Call your voicemail retrieval by dialing the number. Frequently the first speed dial number is used for voicemail. Open the dialer and long press 1 to make the call.

Using a Headset or the Speakerphone

The Speaker button gives you access to the following audio options:



Receiver The receiver speaker is the default speaker you use when you hold the handset to your ear. If another option is in use, return to the Receiver by tapping the receiver icon in the active call screen.



Speakerphone You can direct audio through the speakerphone on the back of the handset. Tap the speaker icon in the active call screen. The speakerphone is designed for one-person use and is not effective as a conference phone.



Headset You can route audio to a headset by plugging a headset into the headset jack on the right side of the handset. The headset jack accepts a four-conductor 3.5mm plug. Use a common headset such as those used by other devices powered by Android; Apple headsets will not work.

Using your PIVOT while you are in a call

You may want to do something else on your PIVOT while in an active call. For instance, you might want to check your calendar or email, open your People app for contact information or browse for an address. You can do any of these activities and continue to carry on a conversation with your active call. To display the active call, tap the home button on the navigation bar and tap the dialer. This opens your call.

Chapter 5: Using the Barcode Reader (8753 only)

Barcode technology enables you to encode and decode information stored in a variety of visual patterns. You can use the barcode reader to input data into any field that accepts text. It is also compatible with the Spectralink Quick Barcode Connector (QBC) program. Contact your system administrator for information about what applications are designed to be used with your barcode enabled PIVOT.

Currently, the Spectralink 8753 handset supports both 1D and a wide range of two-dimensional (2D) barcode types. Two-dimensional barcode symbologies typically encode more data, which may be contained in a smaller barcode label, as compared with 1D barcode symbologies.

1D example (Code 128 symbology)



2D example (QR Code symbology)




Caution: The barcode reader LED

The LED light is not a laser but is very bright. Do not stare at it or shine it into the eyes.

General Operating Instructions

The barcode reader is located at the top rear of the handset. Scanning range for the 8753 is 1 to 18 inches. The reader is ergonomically designed to slant down so that you can comfortably hold the handset in a parallel position without having to twist your wrist to scan a barcode that is on a table.

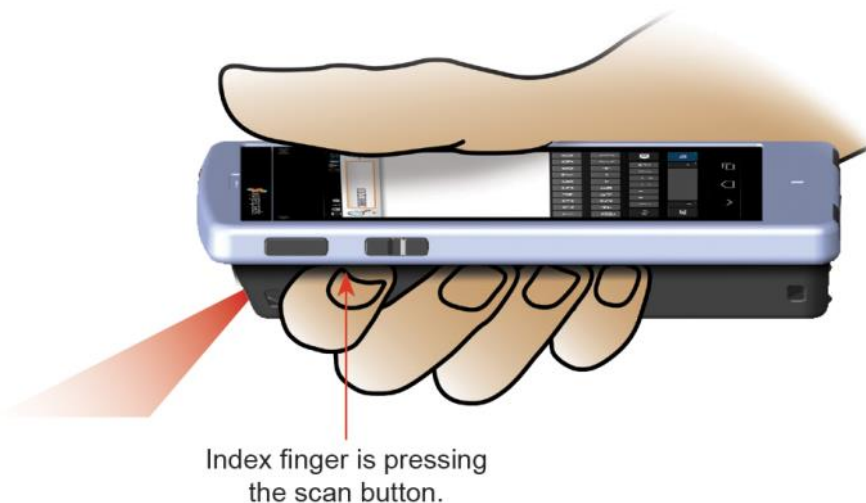

Scanning with the barcode reader:

- 1 Hold the handset so that the barcode reader is pointing to the item you want to scan.
- 2 Press the scan button.

The red and white light comes on. Orient the red line so that it shines across the entire barcode symbol. The lights stay on until you hear a beep (if configured) that indicates the read is successful and then the lights turn off.
- 3 When the red light turns off and a beep or other sound is heard (if configured), the scan is successful and you may release the scan button. If you wish, you can change the barcode reader sound. See [Sound](#).



User Tip: How long does it take to scan a barcode?

The larger the barcode, the more data it contains. 2D barcodes can be very large and contain a lot of information. It can take several seconds to read a very large (e.g. 3000 character) barcode. The light will turn off when the scan is complete and then you may let go of the scan button. If configured, a beep or other sound will also play.

However, if you are scanning rapidly and have a long beep tone, it will not have enough time to play at each scan and will sound intermittent. Set short beep tones for barcode scanning if you anticipate very short time periods between scans and want an audio confirmation. Barcode and Plastic Pipe are examples of short tones. Formalhaut, Ta Da and Procyon are long ones.



Note: Cleaning the scanner window

If you have trouble scanning, clean the scanner window using a monitor wipe or a glass cleaner sprayed onto a soft cloth. Do not spray the scanner window directly. When first using the bar code scanner, ensure that the shipping plastic is removed. The shipping plastic covers the bar code lens and must be removed by its small tab before use, for reliable operation.

Sending Barcode Data to PC Applications using QBC

The Quick Barcode Connector (QBC) enables you to capture barcode patterns and transfer the data to an application running on a host computer. You can think of the system as a wireless barcode scanner connected to a host computer.

Your system administrator will install QBC on one or more computers. Depending on how your administrator sets up the environment, the handsets can transfer data to one computer only (single-endpoint mode), or to any computer in a group of computers (multi-endpoint mode).

In single-endpoint mode, the Spectralink QBC application is installed on one computer. When a handset captures barcode data, the data automatically transfers to the application on that computer.

In multi-endpoint mode, the Spectralink QBC application is installed on a group of computers so that a handset can transfer barcode data to any computer in that group. Although you can connect a handset to any computer in the group, you can connect a handset to only one computer at a time and you can connect only one handset to a computer at a time.

Check with your system administrator to find out which QBC modes are used in your facility.

The QBC application icon in the computer's system tray indicates the connection status. If your handset is not connected to the Spectralink QBC application, the application icon displays a red X. If the handset is connected, the icon displays a green arrow, as shown next.



Barcode scanning in single-endpoint mode

If you are scanning to a PC application in single-endpoint mode, your handset will automatically connect to your computer.

To scan a barcode in single-endpoint mode:

- 1 Open the application on the PC.
- 2 Ensure that your cursor is in the correct position in the application you're using to enter incoming data.
- 3 Scan the barcode(s) you want to enter into the PC application.

Barcode scanning in multi-endpoint mode

To enter barcode data in multi-endpoint mode, you must first connect your handset to a computer.



User Tip: Multi-endpoint status icon

While your Spectralink 8753 handset is set to multi-endpoint mode and connected to a computer, the following status icon is displayed:



To enter barcode data in multi-endpoint mode:

- 1 If you haven't already done so, on the computer, launch Spectralink QBC by choosing Start> Programs> Spectralink> Quick Barcode Connector> Start Quick Barcode Connector.

The application launches quietly and runs in the background. You can confirm that the Spectralink QBC is running by checking for the icon in the status bar of your 8753 handset.
- 2 Move the cursor on the host computer to one of the barcode hotspots, typically a corner of the display.

The application barcode will display temporarily.



User Tip: Displaying the barcode

If you cannot see the barcode display due to interference from other applications running on the host computer, right-click the Spectralink QBC icon in the system tray and choose **Show Barcode** from the pop-up menu. While you have *Show Barcode* selected, Spectralink QBC will persistently display the barcode pattern identifying the host computer until you select **Hide Barcode** from the pop-up menu.

- 3 Use your 8753 handset to scan the barcode on the computer monitor. This connects your handset to the computer.

When your handset has scanned and decoded the barcode pattern and connected to the Spectralink QBC application running on the computer, the following will happen:

- The handset will emit a confirmation tone to notify you that it is connected if it has been configured to make a sound.
- The Spectralink QBC status icon on the handset will appear on the right side of the status icon row.
- The Spectralink QBC status icon on the computer will display in the system tray.
- The computer will display a 'Connected' balloon and a notification sound (if enabled).

The barcode pattern displays the message *Connected* when you move the pointer to the hotspot, as shown next. (If there is no connection between the handset and the host computer, the message *No Connection* displays on the barcode.)



Once you are connected, you are ready to scan data into the computer application.

- 4 Start the application and ensure that your cursor is in the correct position in the application you're using to capture the incoming data.
- 5 Scan the barcode(s) you want to enter into the target computer application.

Disconnecting from the computer

The handset automatically disconnects from the computer when it has been inactive for a while or when a second handset connects to the computer.

Chapter 6: Applications

Two different types of Applications might be available for you to activate. Applications are available from either the Launcher icon or from the Spectralink widget box usually located on home screen 2. Google applications may also be available. Contact your system administrator if you have questions about how to activate and use applications deployed in your facility.



Note: Interactive applications

Certain applications are able to interact with the handset by sending messages or alerts and so forth. Your system administrator will tell you if you have this type of application on your handset. This type of application is not activated by the user.

Spectralink widget box

The Spectralink widget box can be dragged to any home screen. It contains custom applications that have been installed by your system administrator. Simply tap the icon of the application you want to activate. See your administrator if the widget box is empty and you see the message *No configured applications please see your system administrator*.



Note: Don't see your app?

The handset needs to be restarted when the system administrator installs a new app in the widget box. If you do not see an app you expect to see in this box, restart your handset. If it still doesn't appear, contact your system administrator.

All Applications Launcher

Your system administrator may have loaded applications on your handset that are available in the Launcher screen. You can find them in alphabetical order. Just tap the app icon to activate the app. You can also drag an application to any home screen for quick access.



Note: Installing applications

You may not be able to install applications on your PIVOT without specific support from your system administrator. Contact your system administrator if you want to use third-party applications.

Chapter 7: SAFE and Panic Button

Workers can be at risk during security breaches or if personal incidents require immediate attention. The Spectralink SAFE application offers personal monitoring for man down functionality, including “no movement” alarms, “tilt” alarms, “running” alarms and duress calls to an emergency number.



Caution: Features may not be configured

The Panic Button and SAFE applications are available through separate feature licenses and may not be deployed in your facility. Please contact your system administrator for more information.

The Panic Button is the red button on the top of the handset. It provides an alarm and/or instantaneous calling to a pre-programmed emergency number.

The emergency call feature can be programmed to work with both the Panic button and the motion sensors. It can be programmed to activate the speakerphone if an emergency call is placed. Motion applications can be programmed to place an emergency call if the alarm state is triggered.

Panic button and SAFE applications can be coupled with a security alarm application program which can be programmed to receive the alarm and identify the handset, the user and the location of the alarming handset along with other functionality.

Contact your system administrator for full information about how personal safety applications are deployed in your facility.



Warning

The reliability of the Spectralink Panic button and SAFE application depends on the functionality and reliability of the greater infrastructure – the wireless LAN, the LAN, the call server, the central provisioning server, the server hosting location services, the central security system and its servers, the correct configuration of the Spectralink 87-Series handsets, correct installation and configuration management server, and thorough training of personnel.

How the SAFE Application Works

Three conditions of alerting can be activated and each is configurable by your administrator according to the requirements of the facility. The user cannot configure any of these settings. The user may be able to temporarily suspend SAFE motion monitoring if the “snooze” feature is configured.

If any of the SAFE motion conditions occur, the handset will first warn the user about an impending alarm and if the user does not cancel the warning within a configurable number of seconds the handset will start to alarm. These are the SAFE motion conditions:

- “Running” – the handset detects shaking when a user runs for a configurable number of seconds. If the warning is not canceled, an emergency call will be placed, if configured.
- “Tilt” – the handset is not vertical for a configurable number of seconds. If the warning is not canceled, an emergency call will be placed, if configured.
- “No movement” – the handset remains still for a configurable number of seconds, potentially indicating the user is no longer moving. If the warning is not canceled, an emergency call will be placed, if configured.

The SAFE motion conditions are configurable by sensitivity so depending on the sensitivity level configured by the administrator, it may take very little motion/tilt/stillness to trigger a warning or it may take a lot.

If an emergency call is placed, it preempts any active call which is ended by the application.



Caution: Contact your system administrator for training on SAFE motion and Panic button applications

All users need to know what aspects of the SAFE and Panic button applications have been configured and what the result will be if an alarm is activated, either intentionally or unintentionally. The user cannot permanently disable Spectralink Personal Alarms or turn the feature on and off.

Icons and notifications

The Spectralink Personal Alarms icon on the Notification Bar indicates that the Spectralink Personal Alarms functionality is active.

SAFE icons









The shield is used for motion sensors







The check indicates motion sensors are actively monitoring


The closed circle indicates Panic Button is activated.
An open circle indicates Panic Button is not active.

	Charging, motion sensor monitoring paused but Panic Button is enabled and functioning.
	Charging, motion sensor monitoring paused but Panic Button is disabled.

	Charging, motion sensor monitoring disabled but Panic Button is enabled and functioning.
	Motion sensors are monitoring, Panic Button is disabled.
	Motion sensors are monitoring, Panic Button is enabled.
	Motion sensor monitoring disabled but Panic Button is enabled and functioning.
	WARNING: Motion sensor monitoring threshold exceeded. (Panic Button is enabled but not alarming. Note that Panic Button never goes to the warning state.)
	WARNING (with Panic Button disabled)
	ALARMING due to motion sensor trigger.
	ALARMING due to Panic Button press.
	[No icon] Panic Button and motion sensors are not enabled.

When you pull down the Notification Bar you will see more information about the status. For instance:

	SAFE monitoring active, Panic Button enabled Start Snoozing (60)
	SAFE monitoring active, Panic Button disabled Start Snoozing (60)
	Panic Button enabled, SAFE monitoring disabled
	SAFE ALARMING, Panic Button enabled Tilt

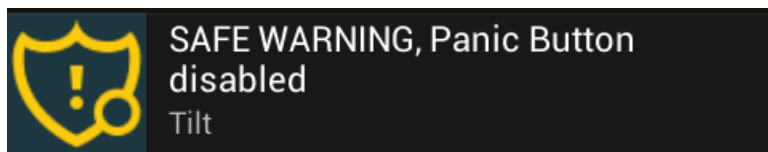



Suspending the motion sensor


If configured, the user can suspend SAFE motion sensing mechanism for a period of time by pulling down the Notification Bar...




...and tapping the notification to begin snoozing:

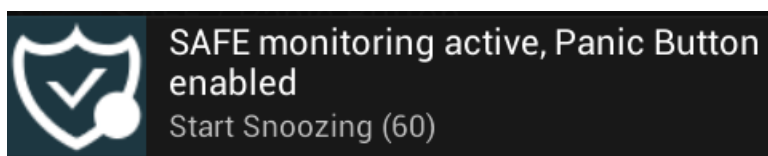

Then the SAFE Notification Bar icon will turn to pause . If you pull down the Notification Bar you will see how much snoozing time is left before monitoring begins again. You can stop snoozing by tapping the notification:

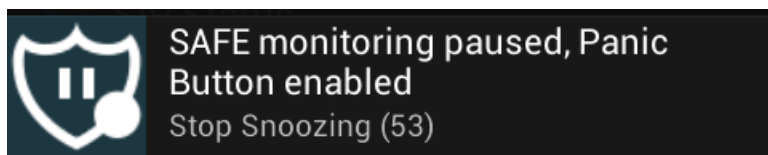

After the Snoozing timeout expires or is manually stopped, the checked shield  appears on the Notification Bar to indicate that monitoring has resumed.




Note: Can't see the snooze option?

If the snooze option is not shown on the notification, use your finger to swipe down on the notification and it should appear.

The Warning state


If a SAFE motion event has been detected, the handset goes into a “warning” state and displays a warning window which shows the type of event(s) and sounds the alert, ramping up the volume as the warning timeout counts down.

SAFE tilt event warning

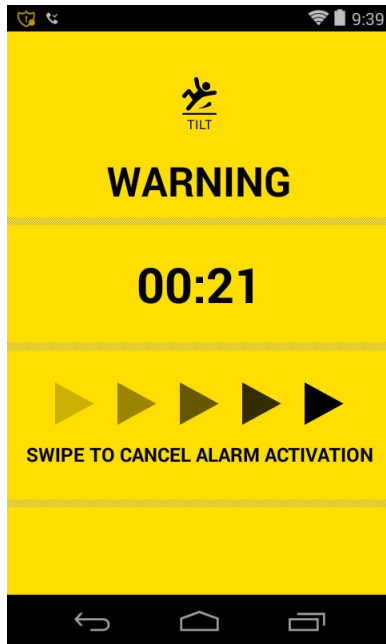


The user can cancel an impending alarm during the warning state by swiping across the screen from left to right over the arrows. The alert audio will turn off and no alarm will be sent to any configured security application and the handset resumes monitoring motion.

While in the warning state, the phone can still be used for any other operation, if desired. However, the warning will continue and can only be cancelled by returning to the warning screen and swiping the arrows. Return to the warning screen by pulling down the Notification Bar and tapping the notification:

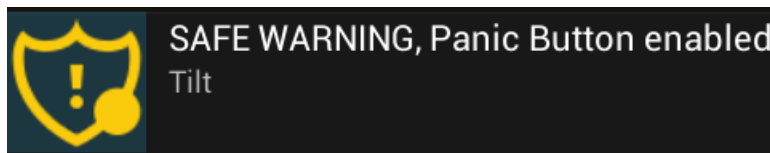

The Alarm state


If the user does not cancel during the SAFE warning state, the handset goes into the alarm state. The alarm sounds at maximum volume and the alarm window appears. If configured, an alarm notification is generated and sent to a security application, and an emergency call is placed.

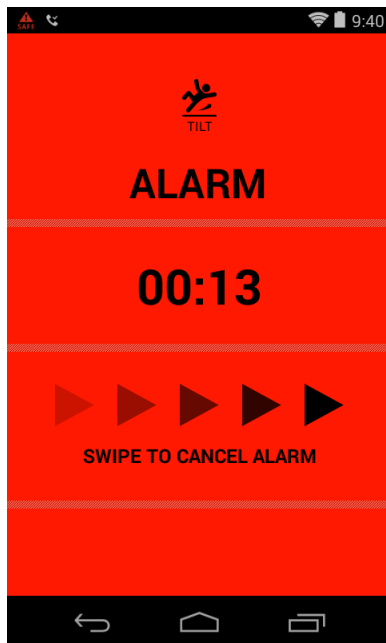

The alarm can be turned down by pressing the volume down button or silenced by using either the volume down button or the power button. The alarm state will continue until cleared.

Once an alarm state is reached, the user cannot retract the alarm but can cancel the alarming state and return the handset to the monitoring state by swiping across the screen from left to right over the arrows. Cancelling the alarm returns the handset to the monitoring state. If an emergency call has been placed because of the alarm then the call must be terminated, as would any normal call.

While in the alarm state, the phone can still be used for any other operation, if desired. However, the alarm will continue and can only be cancelled by returning to the alarm screen and swiping the arrows. Return to the alarm screen by pulling down the Notification Bar and tapping the notification:

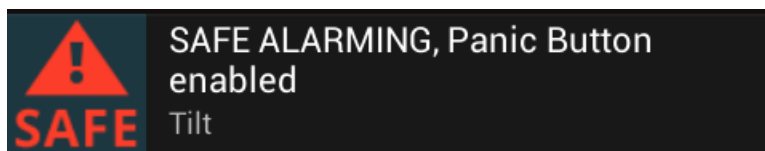

Using the Panic Button

When you press the red emergency button on the top of the handset, the alarm state is entered. A loud warning tone may sound and a call may be placed to an emergency number.

Panic Button icon


The Panic Button icon is indicated by a closed circle on the SAFE shield. An open circle indicates that the Panic Button is not configured.

Activating the Panic Button

The Panic Button can be configured in several different ways. Contact your system administrator to find out which of these setting is active in your application:


- Is the Button activated by a long press, two short presses or either?
- Does pressing the Panic Button place an emergency call?
- Does pressing the Panic Button sound an alarm through your handset's speakerphone?
- When an emergency call is placed, does the audio come through the speakerphone?

If you press the Panic button in a manner that is not consistent with the configured trigger, a popup toast will tell you how to successfully trigger a panic alarm.

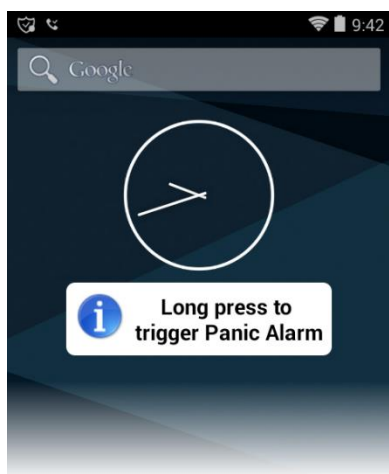

An emergency call preempts and terminates any existing active call. The emergency call is ended by hanging up the call, as with a normal call, or by the far end hanging up. However, the handset remains in the alarm state until manually cleared by the user by swiping across the screen from left to right over the arrows.

Panic button alarm state

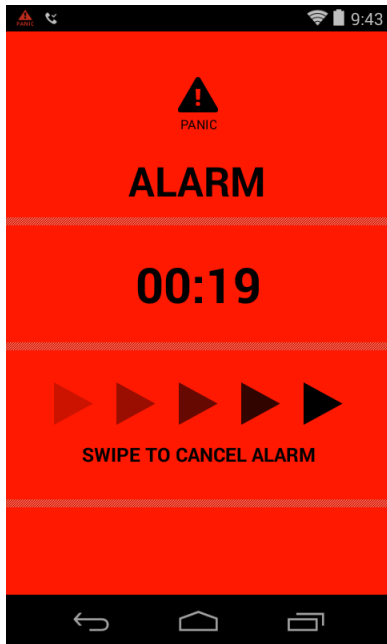


Chapter 8: Push-to-Talk (PTT)

Spectralink 87-Series handsets can operate in a group broadcast mode called Push-to-Talk (PTT). With PTT, the administrator subscribes users to some or all PTT channels to receive, and optionally send, PTT broadcasts.

What is PTT?

In PTT mode, the handsets behave like walkie-talkies: your handset transmits audio over a built-in speakerphone and recipients can respond to your message. Your system administrator defines which PTT channels are subscribed to by your handset. If the system administrator has disabled PTT you will get a disabled notice when you launch the PTT application.

PTT activity is shown on the status bar and you can open the app to see which channel is active. For the fastest access to PTT the app, drag the PTT app icon to your Home screen.

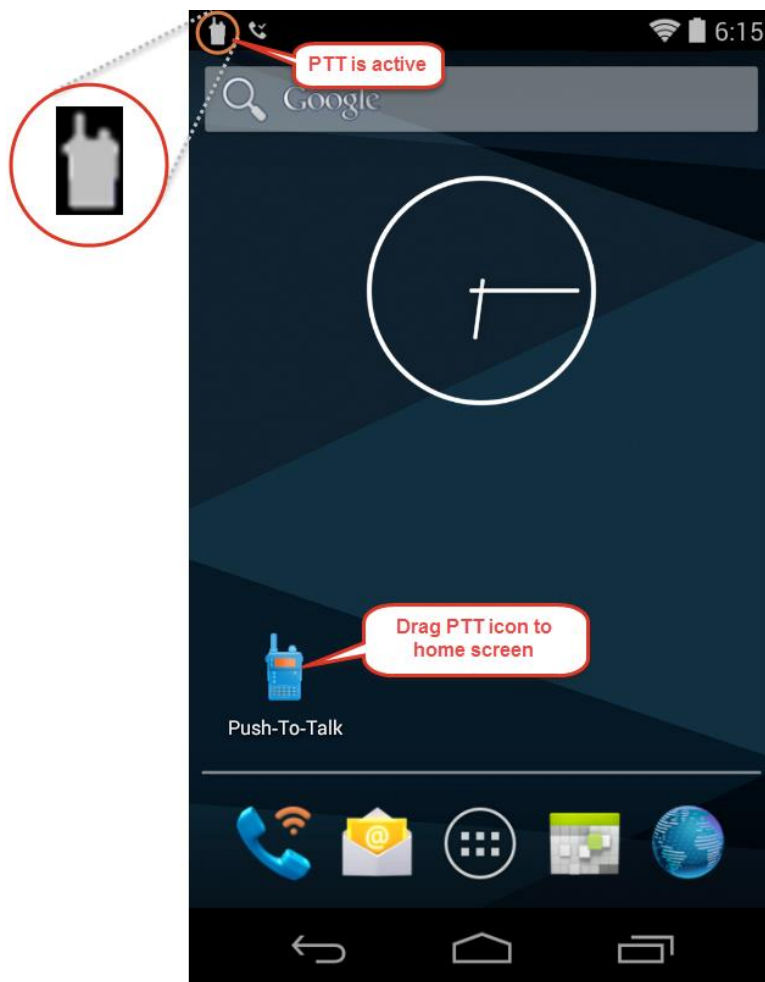

There are three PTT channel types:

- **Normal Channels**

PTT channels 1 through 23 can be used as normal channels by any subscribed handset user. The administrator can create a label, comprised of up to 15 characters, for channels 1 to 23 to help users identify PTT channels. The default label for channel 1 is ALL.

Normal channels are usually allocated to various groups of users and recipients. Examples include channel 2 for the maintenance group, channel 3 for Hospital Admitting, channel 4 for the IT Helpdesk.

- **Priority Channel**

When you are subscribed to the priority channel (channel 24) you will hear transmissions on the priority channel unless you are in a phone call or an emergency PTT conversation is already playing. (i.e. Priority channel PTT conversations interrupt normal channel PTT conversations.) The priority channel plays transmissions at the currently set audio level.

- **Emergency Channel**

When you are subscribed to the emergency channel (channel 25), your handset will hear emergency transmissions unless you are on a phone call in which case you will hear an in-ear alert and the screen will display an option to hold the call and listen to the broadcast. Emergency transmissions play out on the speaker at an administrator configurable audio level. The audio level will return to normal for subsequent non-emergency transmissions.

Subscribing and Unsubscribing to Channels

Contact your administrator to change your subscriptions and transmission permissions. Only the channels that you are subscribed to display in your PTT app. Only the system administrator can change what channels display.

Using the Push-to-Talk Feature

The PTT Application




Note: Does my handset need to be unlocked?

Some facilities require that handsets be locked when they are idle. You can receive PTT broadcasts on your locked handset but in order to respond or initiate your own PTT broadcast, you must unlock it.

To open the PTT app

- Press and quickly release the **Talk** button
- Tap the PTT application icon in the Launcher or on one of your Home screens.



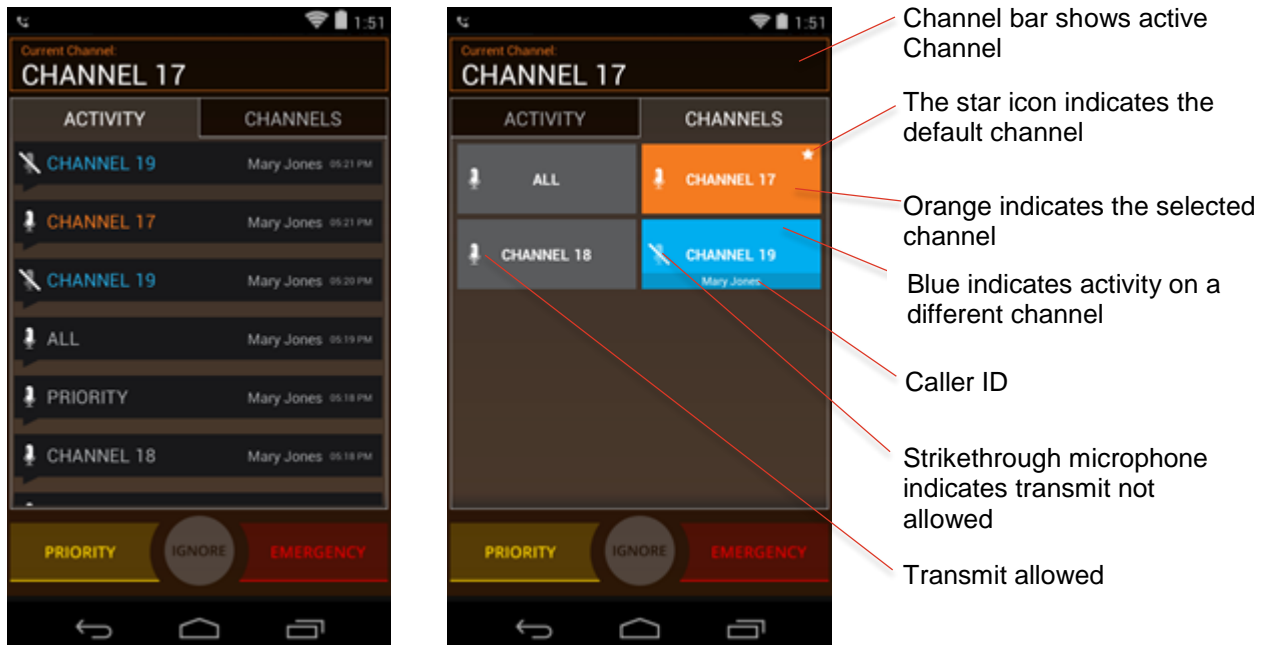

The PTT app displays the current channel and either the Activity tab or the Channels tab, depending on which tab was last viewed.



Key:

Color	Activity tab	Channel tab	Description
Orange	text	background	Indicates the currently selected channel (as shown in the current channel bar). When you press the Talk button, this is the channel you will transmit on.
Blue	text	background	Indicates a channel other than your selected channel that has current PTT activity. You cannot hear or contribute to this activity unless you switch to this channel.
Gray	text	background	Indicates a channel other than your selected channel that no longer has current PTT activity (Activity tab) or is currently inactive (Channels tab).

Icon	Description
	Default channel icon
	Can Transmit icon - You are allowed to transmit on this channel.
	Listen Only icon - You are not allowed to transmit on this channel

<i>Icon</i>	<i>Description</i>
	Receiving icon – you are receiving a transmission on this channel
<i>Button</i>	<i>Description</i>
	Priority Channel Selector
	Ignore Button
	Emergency Channel Selector

What is a PTT broadcast?

PTT broadcasts are similar to chat room sessions. You can either start a new conversation on an idle channel or join a conversation on an already active channel. Your default PTT channel is probably the channel that you use for most conversations.

PTT uses the term Current Channel to refer to the channel that you are currently “tuned” to – meaning it is the channel you are listening to conversations on and it is the channel you will transmit on if you press the Talk button. The PTT app makes it very easy to select the channel you wish to talk on – which will become the Current Channel.

Transmitting a PTT broadcast

To start a PTT conversation

To transmit on the current channel, the handset must be unlocked.

- 1 Press and hold the Talk button.
It takes about 2 seconds for the transmission to begin on the other end which will hear a beep to indicate that a broadcast is incoming.
- 2 Hold the handset microphone approximately two inches from your mouth and speak while holding the button.
You will notice the microphone icon show up in the Channel Bar when you are transmitting. Do not let go of the button until you are finished talking.
- 3 When you let go of the button, the channel enters a ten-second wait period during which the channel is open to receive or send additional broadcasts.

If no activity occurs within the wait period, PTT session is closed. The wait period keeps the channel open and gives any handset subscribed to that channel time to join the conversation.

To transmit on a different PTT channel

1 If there is an active PTT session on a different channel than the channel you want to talk on, open the PTT app and click the channel from the Channels tab or click a leg of a conversation on the channel from the Activity tab. The newly selected channel will show in the Channel Bar.

2 Press and hold the Talk button. Hold the handset microphone approximately two inches from your mouth and speak.

You will notice the microphone icon show up in the Channel Bar when you are transmitting.


3 When you release the Talk button, the handset enters a ten-second wait period during which the channel is open to receive or send additional broadcasts.

If no activity occurs within the wait period, the channel is closed. The wait period keeps the channel open and gives any handset subscribed to that channel time to join the conversation.

To transmit on the priority channel:

1 Tap the priority channel selector.

2 Press and hold **Talk**, and deliver your PTT message.

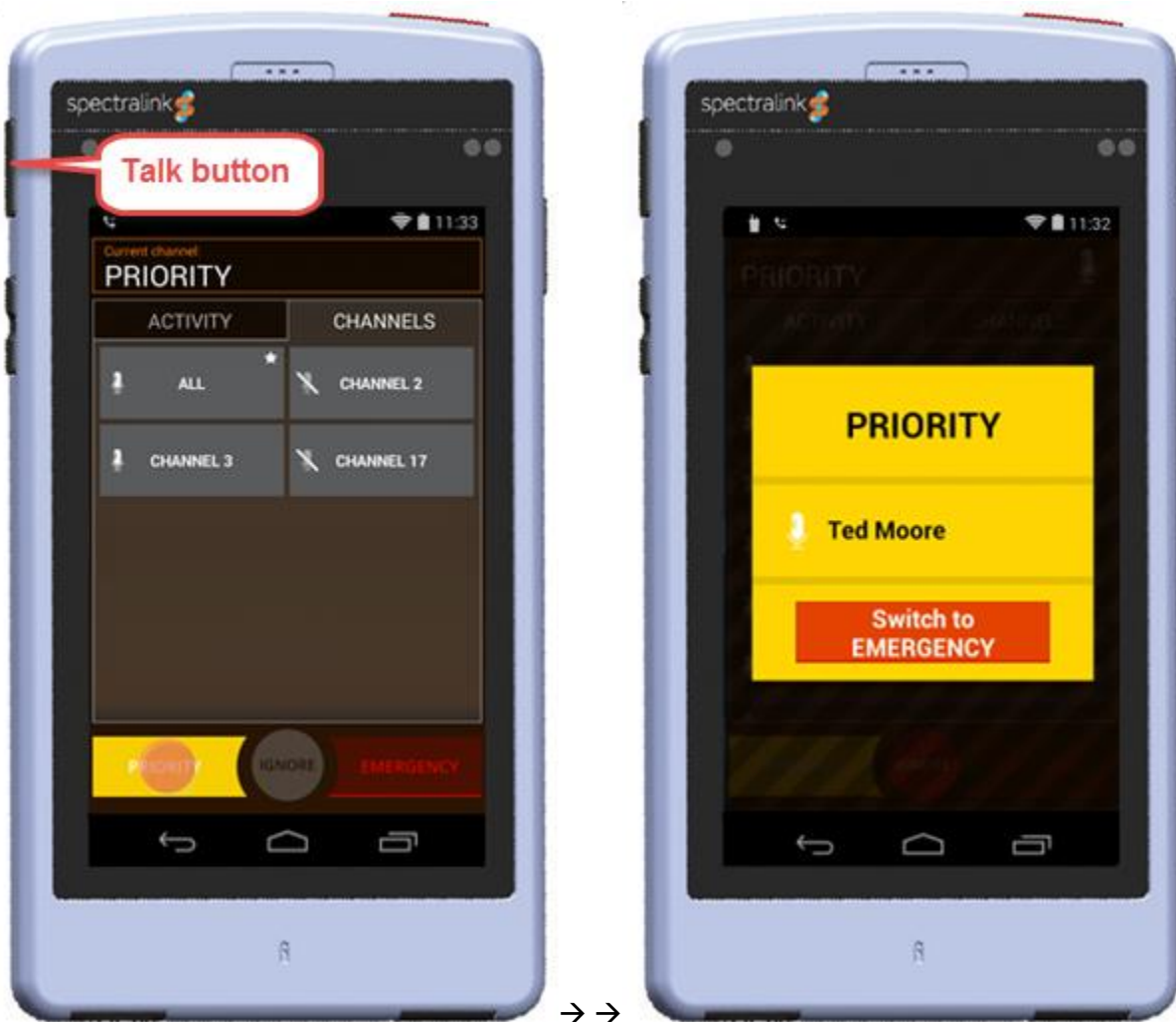

Your transmission will interrupt all other PTT activity. If the recipient is in a non-emergency phone call, an in-ear alert will play to notify the user of the priority broadcast. You may tap the Switch to Emergency button to broadcast on the Emergency channel instead.

To transmit on the emergency channel

- 1 Tap the emergency channel selector.
- 2 Press and hold **Talk**, and deliver your PTT message.



Your transmission will interrupt all other PTT activity. If the recipient is in a non-emergency phone call, an in-ear alert will play to notify the user of the emergency broadcast.

Receiving PTT Broadcasts

Your handset will receive PTT transmissions differently depending on the priority of the PTT Channel whether you are on an active call.

PTT channel priority

The handset will alert you to PTT transmissions from normal channels, the priority channel, and the emergency channel.

Normal channel transmissions

- The PTT icon will appear on the status bar.
- If you are not in a call, you will hear the PTT beep and the audio.

Priority channel transmissions

- The yellow Priority display will appear on the screen.
- If you are not in a call, you will hear the PTT beep and audio.
- If you are in a call, you will hear an in-ear PTT alert but not the PTT audio unless you open the channel using the notification.

Emergency channel transmissions



- The red Emergency display will appear on the screen.
- If you are not in a call, you will hear the PTT beep and the audio will be at max volume.
- If you are in a call, you will hear the PTT beep at max volume and the display will give you a Hold Calls button that you can tap to hear the emergency broadcast.
- If you are on an emergency call, you will hear an in-ear PTT alert.

To view PTT activity

If you have enabled notifications, you will hear a beep to indicate an incoming PTT broadcast and the PTT icon will display in the status bar:



- You can pull down the notification drawer to find out about the PTT activity. Tap the notification to open the app. The example below shows an incoming call on channel 17.

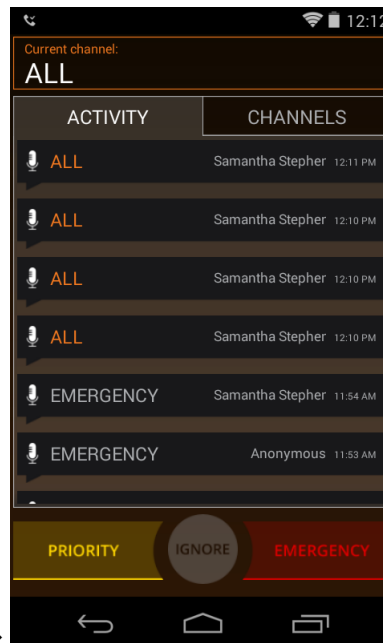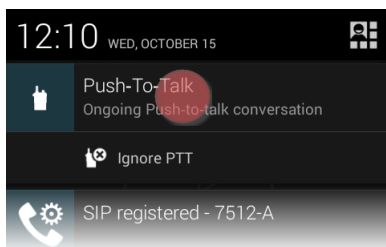

To answer a PTT broadcast


- Press and hold the **Talk** button while the channel is idle. Only one person can broadcast at a time.

To change the PTT volume


- Press the Volume up or Volume down buttons on the side of the handset while in an active PTT call.

To ignore a PTT broadcast

- Pull down the navigation drawer and tap the Ignore option to prevent PTT interruptions for the period that this particular broadcast remains open.

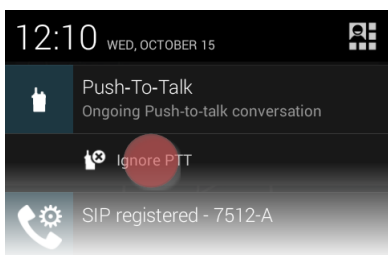

- You can also open the app and tap the Ignore button.

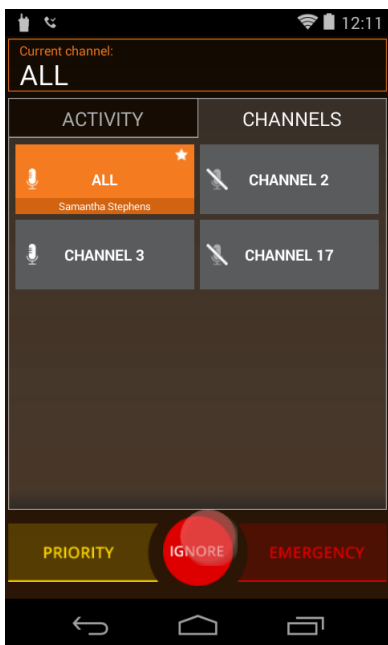

Receiving PTT Broadcasts while in Call

The handset will receive PTT broadcasts differently depending on whether or not you are in an active call. If you are in a phone call you will not be interrupted by a routine PTT call. No beep will sound.

To listen to a PTT broadcast

- Swipe down the notification drawer and tap the notification. This joins you to the PTT broadcast.

To answer a PTT broadcast

- Swipe down the notification drawer and tap the notification. This joins you to the PTT broadcast.
- Press and hold the Talk button while the channel is idle.

To change the PTT volume

- Press the Volume up or Volume down buttons on the side of the handset while in an active PTT call.


To answer a telephone call during a PTT broadcast

- When the incoming call screen is presented on the display answer the call in the normal manner. The PTT audio will be muted.

To start a telephone call during a PTT broadcast

- Switch to the Phone app and start the call in the normal manner. The PTT audio will be muted.

Receiving Priority PTT Broadcasts

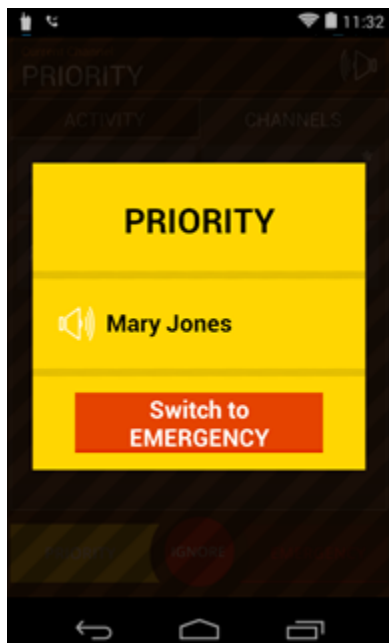

Receiving Emergency PTT Broadcasts

When idle



When in call



The volume slider appears when you receive an emergency call when the handset is idle. You can use the volume slider to lower the volume, if desired.


When you are in a telephone call, you will hear an in-ear alert. Use the Hold Calls button to put your active call on hold and hear the Emergency broadcast.

Changing the Default PTT Channel

The default PTT channel is the one that is transmitted on when you press and hold the Talk button if there are no ongoing conversations on another PTT channel. You can change your default channel on the Channels tab on the PTT UI.

To change the default PTT channel

- 1 Open the PTT app and click on the Channels tab.
- 2 Tap and hold the channel button you want to be your default channel.
- 3 Tap the popup “Set as Default Channel”.
- 4 You will see the default channel star icon on the selected channel and it will display an orange background (because it is selected).

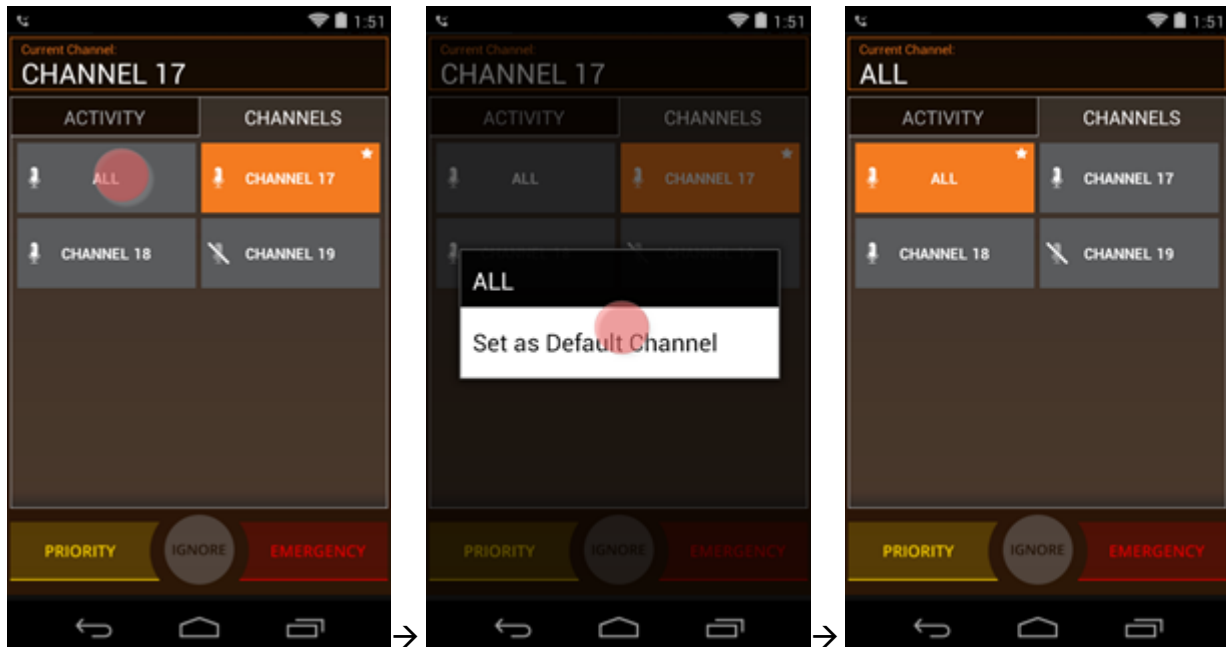

Chapter 9: Chargers and Battery Packs

Your handset is powered by a removable Battery Pack that you will need to fully charge before first using the handset and then recharge periodically. You may charge the Battery Pack while it is removed from or while attached to your handset. Your system administrator will probably supply a USB charger with your handset.



Caution: Battery compatibility

Only use Spectralink 87-Series Battery Packs with Spectralink 87-Series Handsets. Check the label on the Battery Pack and ensure the model number is DM351.



Note: Battery life

How intensively you use your PIVOT determines how frequently the Battery Pack will need to be charged. The handset is designed for 20 hours of standby time. Under moderate use circumstances, a Battery Pack will last a full 8-hour shift.

A good rule of thumb is to plug your handset into your USB charger after every shift or overnight. Remember to leave it on. A handset that is powered down will not charge.

Under intensive use, Battery Pack replacement may be needed during a normal shift. Facilities frequently provide a charging bank with extra Battery Packs for replacing Battery Packs during the day if intensive use is usual.

Contact your system administrator if you need to change Battery Packs during your shift.

How usage affects Battery Pack life

Certain handset features require more battery capacity than others. Any battery life estimates are highly dependent on phone usage and the features that are deployed as well as the environment itself. Battery life varies according to the battery power that is required to support high draw features. Having a removable Battery Pack allows you to switch to a fresh Battery Pack whenever you hear the low battery signal.

PTT

When Push-to-talk is deployed, Battery Pack requirements are significantly increased depending on the amount of actual broadcast activity. How significant an increase depends on how often any given handset receives a broadcast, how long the broadcast channel remains open, and how often one initiates a broadcast.

Barcode scanning

Some facilities make frequent use of the barcode scanning feature. When you press the scan button, it activates the scan LED which puts a high demand on the Battery Pack. Stay aware of this and use your barcode scanning efficiently.

Background brightness and wallpaper

Like any device with a display, the brightness of the display, the background and the length of time the display is active are all factors that contribute to battery life. If you set an active background, make the display very bright and keep it active for the maximum time possible, your Battery Pack will have shorter times between charges than if you use more conservative settings. Adjust your Background settings by going to **Settings> Display** and selecting the area you want to adjust.

Excessive roaming

Your handset is designed to roam efficiently throughout a well-planned facility. However, some places have “dead spots” or areas where the access points do not provide adequate coverage. In these weak signal areas, the handset is continuously roaming trying to find a better connection so you can make and receive calls. If your facility has these areas, alert your system administrator. If you are out of range or on the edge of coverage, turn your handset off to preserve battery life.

Feature combinations

The battery life factors are cumulative. If you are using several features, you will probably need to replace your Battery Pack throughout the day.

Charging the Battery Pack

You must fully charge your handset’s Battery Pack before you use it for the first time. Depending on the type of charger you use, this might take several hours.



Caution: Battery compatibility

Only use Spectralink 87-Series Battery Packs with Spectralink 87-Series Handsets and chargers. Check the label on the Battery Pack and ensure the model number is DM351.



Warning: Shock risk

Risk of electric shock. Do not expose charger or battery to liquid, vapor, or rain.

The Battery Pack icon at the upper right corner of the Status Bar of your handset is completely full when it is fully charged. When the handset is in a charging device and charging, the Battery Pack icon has a lightning bolt symbol through it.



Fully Charged Battery Pack



Battery Pack connected to power source

Charger Options

Two types of Chargers are available. Contact your system administrator if you need assistance with either of these options. The USB charger charges a single handset through the USB port on the bottom left side of the handset. The Quad Charger can charge up to four Battery Packs at a time. The Battery Packs for the Spectralink 87-Series Wireless Telephone can be charged with either Spectralink 87-Series charging options.

Using the USB Charger

The USB Charger (UE08WCP-050150SPA), shown next, plugs into the USB port on the bottom left of the handset. Plug the charger into a wall socket. Carefully open the hinged rubber cover on the bottom left side of the handset and plug in the USB cable.



Using the Quad Charger

The Quad Charger is used in facilities where the USB charger is not efficient and extra Battery Packs are needed during a shift. It can charge four Battery Packs simultaneously and can be mounted on a horizontal or vertical surface. Three Quad Chargers can be used in series to form a charging bank. See *Spectralink 87-Series Quad Charger Quick Start Guide* for mounting instructions.

The Battery Pack slides into the charging slot as shown in the photo below. The LED beside the occupied bay indicates the charging status of the Battery Pack in the bay. A discharged Battery Pack will take approximately 3 hours to fully charge.



The following table lists the LED indications and charger status.

Quad Charger LED Indicators.

LED Indication	Battery Pack charge LED
Solid red	Charging
Solid green	Fully charged
Off	Not properly seated or empty
Blinking red	High temperature or error

If there is an error and the LED is blinking red, the LED will continue to blink red—even when you remove the Battery Pack—until a Battery Pack without error is placed in the slot. Power cycling the unit will also clear the error indication



Caution: Overheated Battery Packs

If an overheated Battery Pack is detected, charging will suspend and the LED will blink red. Remove the Battery Pack and let it cool before attempting to charge it again.

About Battery Packs

The lithium ion rechargeable Battery Pack gives you approximately 20 hours of standby time. Standby time occurs when the handset is on but not in an active call. Standby and charging times vary according to the type of charger you are using.

A discharged Battery Pack will take approximately 3 hours to fully charge.



User Tip: Fastest charging

The Quad Charger is the fastest way to charge a Battery Pack. If you use the handset while it is charging with the USB Charger, charging time will be prolonged.



User Tip: How to determine type of battery

The Battery Pack type is printed on the label. Remove the Battery Pack and view the label to determine the type of Battery Pack you have been issued.



Warning: Explosion risk

WARNING: THERE IS A RISK OF EXPLOSION IF YOU REPLACE THE BATTERY PACK WITH AN INCORRECT TYPE. ONLY USE SPECTRALINK MODEL DM351 WITH SPECTRALINK 87-SERIES WIRELESS HANDSETS. DISPOSE OF USED BATTERIES BY RETURNING THEM TO THE DISTRIBUTOR OR PRODUCER WHO SUPPLIED THEM.

Indications of Low Battery

The Battery Pack icon in the status row on the handset's screen indicates the amount of charge remaining.


Battery Pack Removal and Replacement

PIVOT 8741/8753


To remove the Battery Pack from the handset, slide the carrier latch towards the phone body. The Battery Pack releases outward.



To replace the Battery Pack, slip the tabs at the top of the carrier into the corresponding slots at the top of the battery cavity.

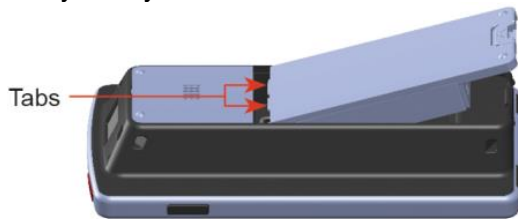

Gently press the bottom of the Battery Pack until it snaps into place. The Battery Pack should snap into position with minimal pressure.

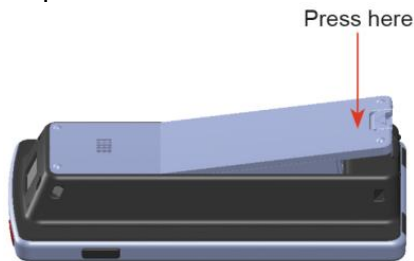

PIVOT:S 8742

To remove the Battery Pack from the handset, slide the carrier latch towards the phone body. The Battery Pack releases outward.



To replace the Battery Pack, slip the tab at the bottom of the carrier into the corresponding slot at the bottom of the battery cavity.



Gently press the top of the Battery Pack until it snaps into place. The Battery Pack should snap into position with minimal pressure.



Battery Safety Notices



Warning: Proper disposal of Battery Packs

Do not throw away the Battery Pack with your domestic waste. Return used Battery Packs to your supplier or servicing agent.



Warning: Handle batteries with care

All batteries can cause property damage and/or bodily injury, such as burns, if a conductive material such as jewelry, keys, or beaded chains touches exposed terminals. The conductive material may complete an electrical circuit (short circuit) and become quite hot. Exercise care in handling any charged battery, particularly when placing it inside a pocket, purse, or other container with metal objects.



Warning: Risk of fire or explosion

WARNING— The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not disassemble, heat above 45°C, or incinerate. Replace battery with Spectralink model DM351 only. Use of another battery may present a risk of fire or explosion.

Notes on Battery Packs

- Battery Pack life will vary depending on your handset's model and its features and system infrastructure.
- Ensure that you fully charge the Battery Pack the first few times, as the Battery Pack becomes most efficient after the first few charge cycles.
- If multiple Battery Packs are supplied with your handset, Spectralink recommends that each be fully charged upon receipt to prolong battery life. Battery Packs will slowly lose charge if unused. To maintain battery potential, charge unused Battery Packs occasionally or alternate Battery Pack use.
- After a length of time, Battery Packs will lose the ability to maintain a charge and to perform at maximum capacity and will need to be replaced. This is normal for all batteries.
- If your Battery Pack requires replacement, contact your system administrator.
- Charging of a battery in the handset can only be done while the handset is turned ON. If the handset is turned OFF it will not charge.
- If the handset does not charge, clean the Battery Pack, the charger, and the handset contacts with an alcohol swab.
- Any battery that exhibits swelling, cracking, or other abnormality should be disposed of promptly and properly. Do not use a battery in this condition.

- Dispose of any used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire.
- The handset is not designed for and will not function properly with third party batteries.

Chapter 10: Sync Calendar and Email

Your system administrator may have enabled your handset to sync with the Exchange Server in your facility that hosts your Outlook calendar and email. If so, you can set up the sync by following these steps:




Caution: Steps may vary!


Your system may operate differently from the steps outlined here. Contact your system administrator for exact instructions for syncing your email and calendar.

To sync your calendar and email


- 1 Tap the Email application icon.
- 2 Enter email address and password and tap Manual setup.

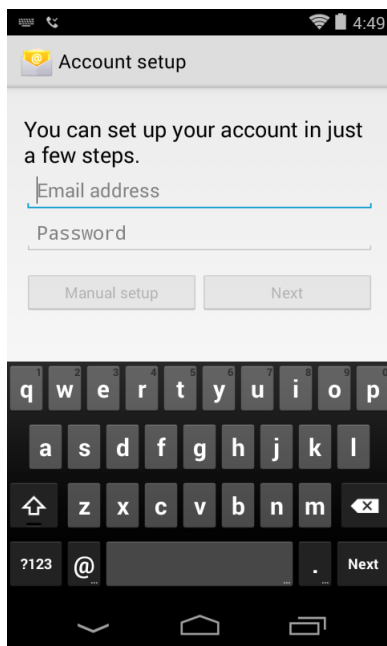

- 3 Tap Exchange.

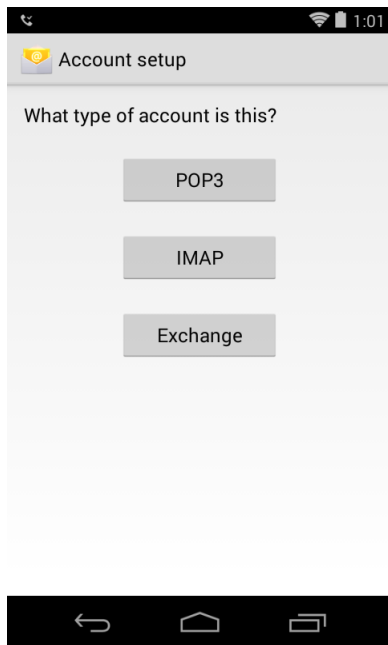

- 4 Modify the Domain\Username field to show your username as directed by your system administrator. Usually this is the email address as entered in step 1. The password should be populated with the one already entered in step 1.

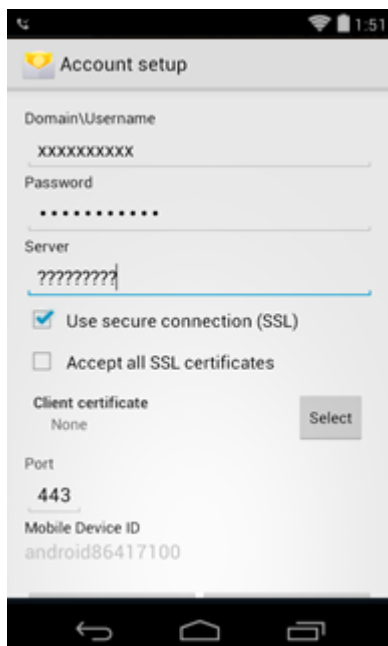


- 5 Enter the server name as directed by your system administrator in the Server field
- 6 Check "Use secure connection (SSL)" as directed by your system administrator.
- 7 Tap Next.




Caution: Consult with your system administrator

Complete the next steps only after receiving specific instructions from your system administrator.


- 8 Tap “OK” on Remote security Administration warning.

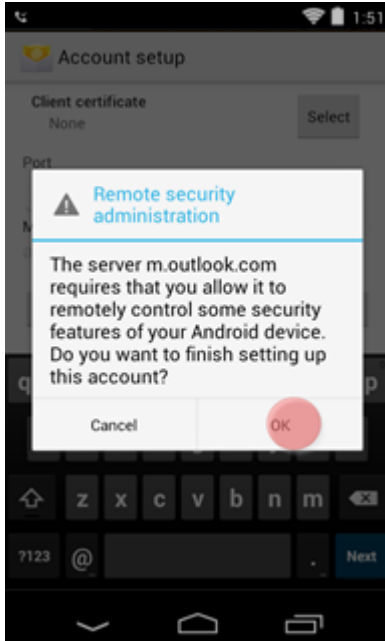


- 9 Leave the settings on the screen as default and tap “Next”

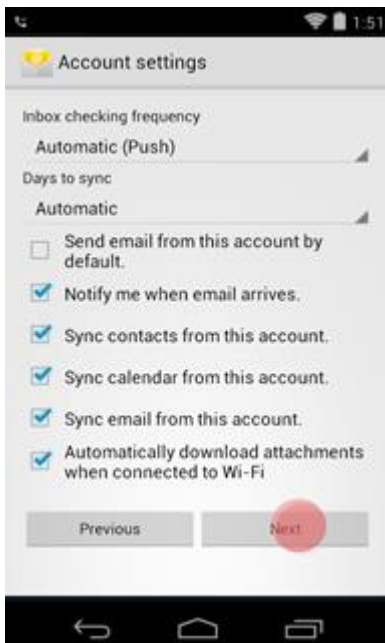


- 10 Change account name if you wish and tap “Next”.

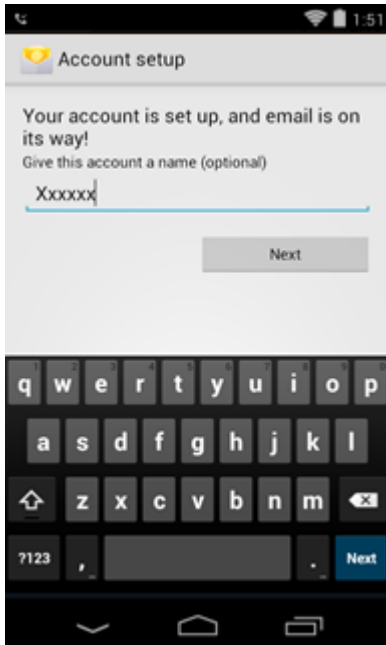


11 Tap “Activate” on Active device Administrator.

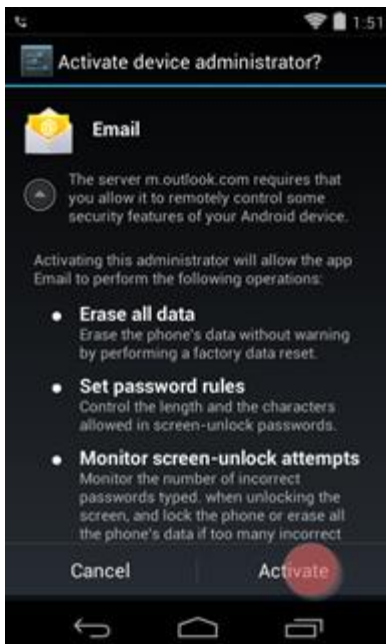

12 Usually screen security must be set. Tap PIN or Password, enter your code and tap Continue. You will be prompted to confirm the code.

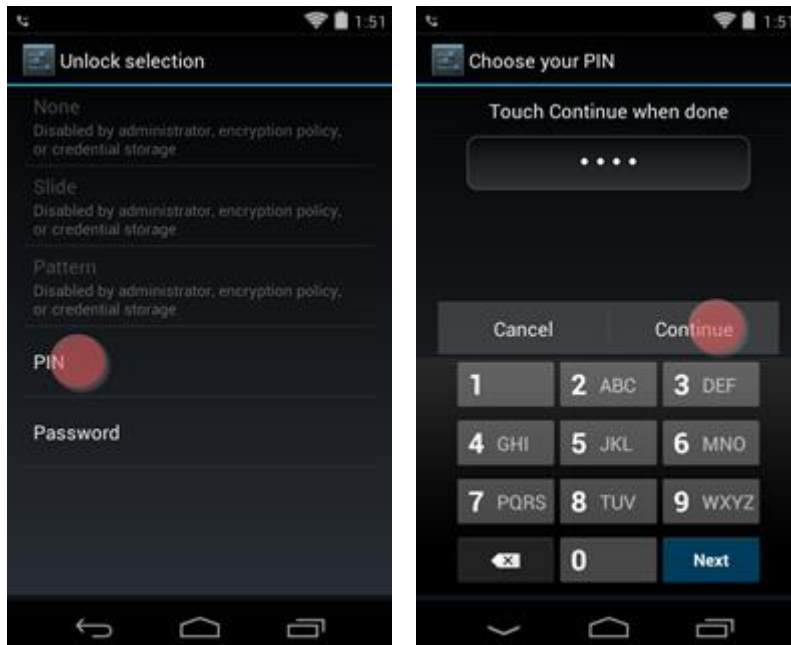

- 13 You should see “Waiting for sync” screen and you may get a number of notification alerts.

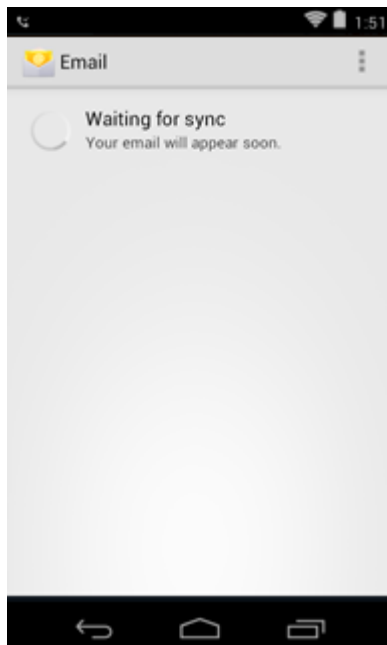


- 14 Once the sync is complete, you will see the lock screen with the email icon in the status bar if you have new emails. Enter your PIN or password to open your email inbox.
- 15 Now you can open your inbox from the notification drawer (if you have new emails) or from the Email icon. Your Calendar and People app will also be populated from your data in Outlook.

Chapter 11: Settings

Use the Settings application to manage things like ringtones, language, etc. and to find information about your PIVOT.

Settings Options

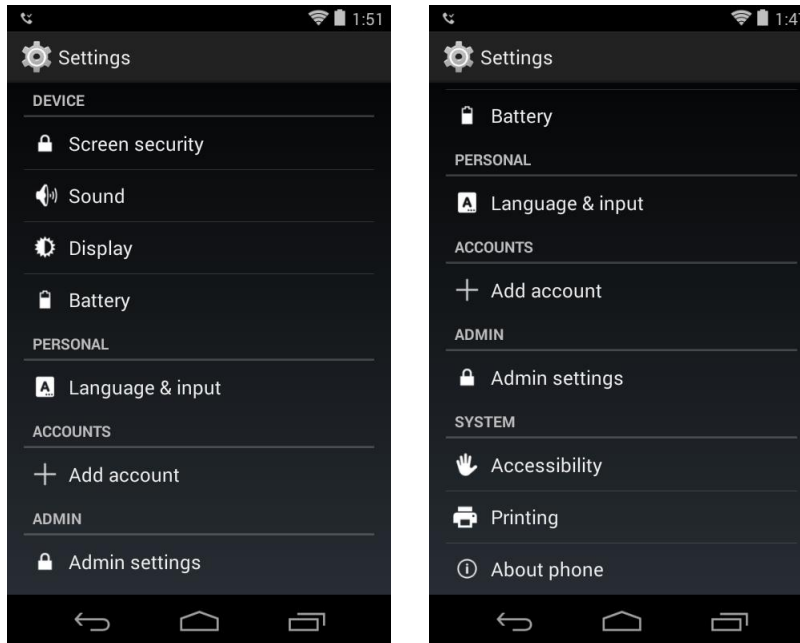


Sound

Change volumes, tones, and other sound settings:



Note: Web API volume setting

The Web API volume setting **Settings > Sound > Volumes > Web API volume** controls the loudness of any alert that you might receive if your handset connects with call systems or the like.

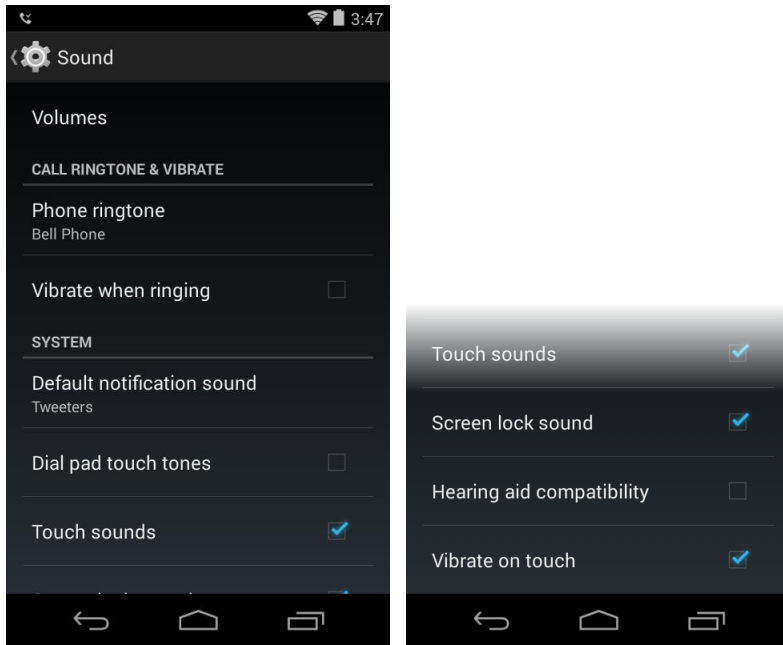


Display

Adjust screen brightness, change wallpaper and other options.



Battery

Check the status of the battery, systems that are using battery life and watch it charge and discharge on this screen.

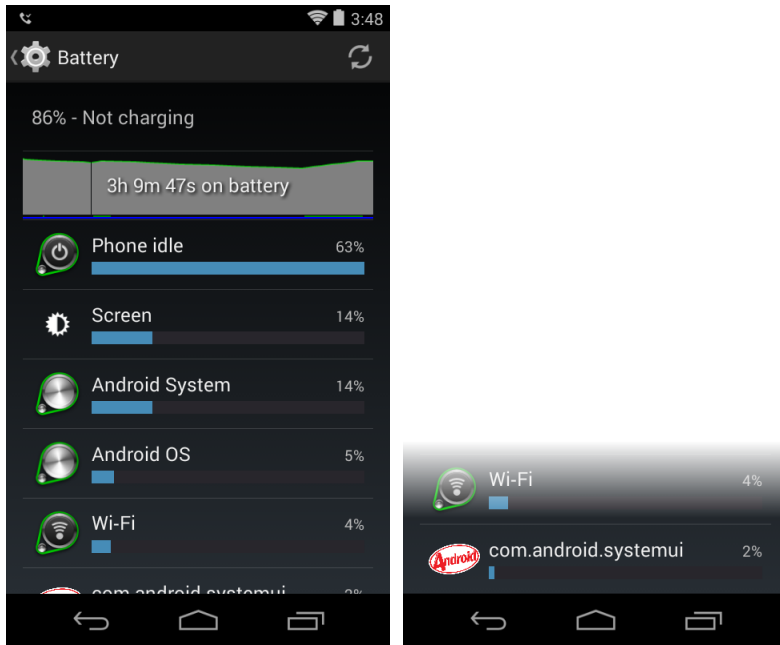


Language & input

Change the language and administer other input tasks.




Note: Changing the language

The Spectralink phone application is customized for your PIVOT and combines standard Android platform functionality with unique Spectralink functionality. The unique-to-Spectralink text is translated into French, German, Italian and Spanish. If you chose a different language, the unique text will not be translated and will appear in English only.

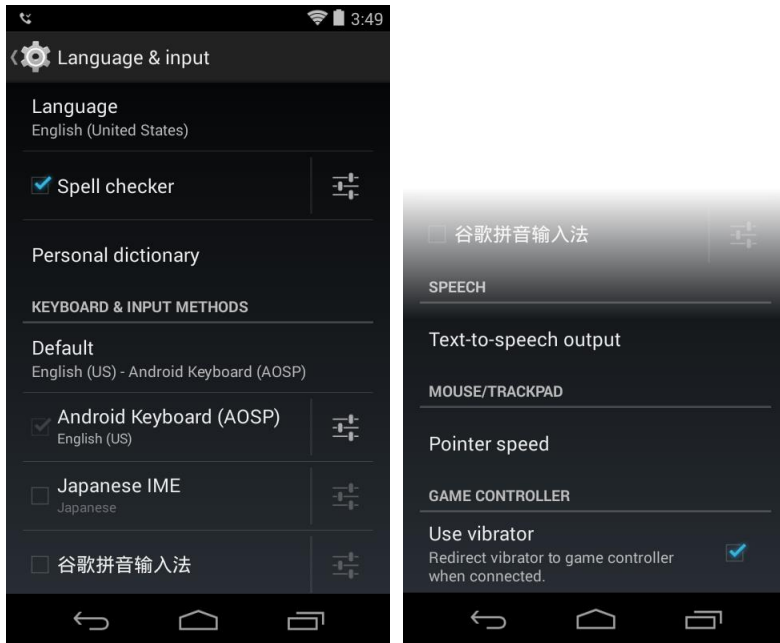


Add an account

Allows you to set up different accounts.

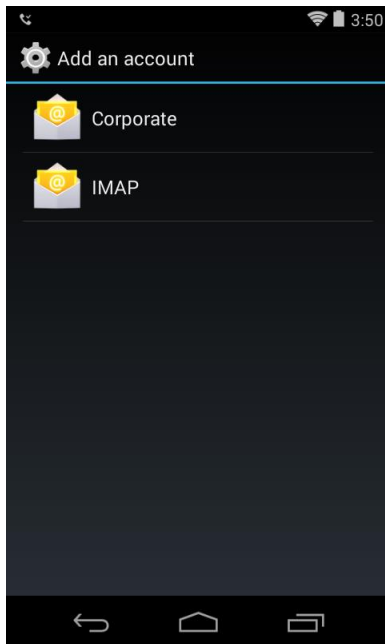


Accessibility

Accessibility options allow you to customize the interface to accommodate special needs.

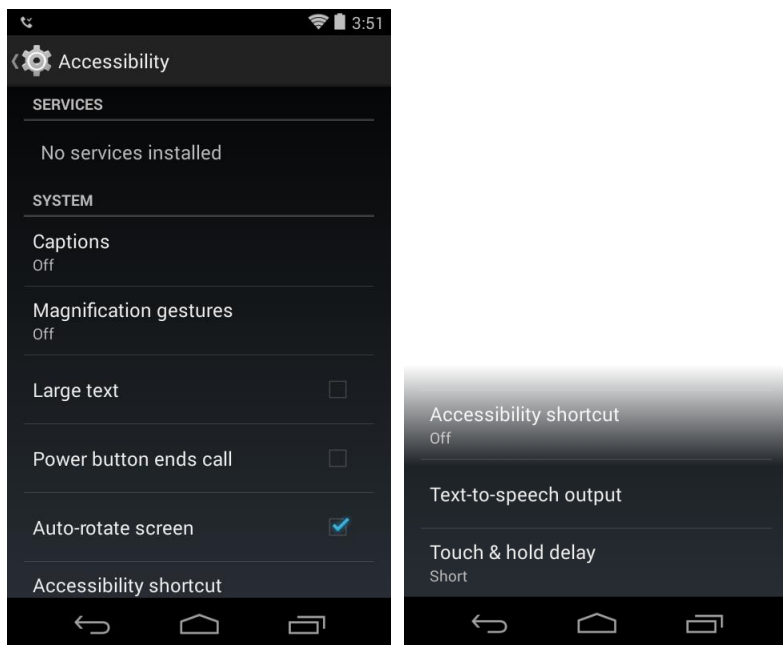

About phone

The About phone screen gives you technical information about the handset.

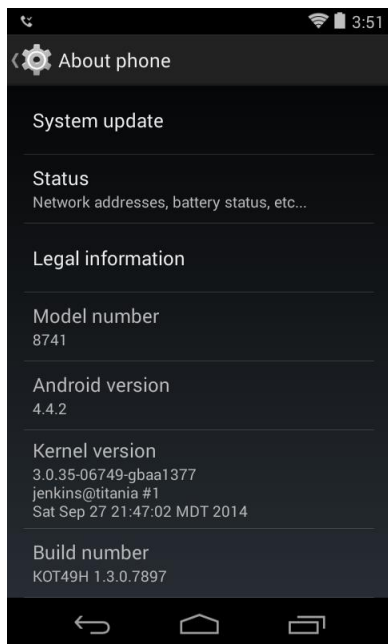


Updating Handset Software

From time to time, your system administrator may provide new or updated software for your PIVOT. If the handset has found new software to install, the download icon will display in the status bar.



Update the software

- 1 When the download is complete, you can swipe-down the notification drawer which will display a notification that the update is ready.

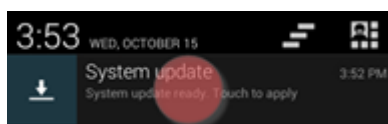

- 2 To apply the update: Tap the System update notification to apply the update and reboot the handset.

Alternately, the handset will apply an update when it is power cycled or when the Battery pack is replaced.

Verifying the update

The system administrator may ask you to check to see if the software successfully installed. Navigate to **Settings> About phone> Build number** and note that the build number matches the number given by the system administrator.

Chapter 12: Other Things You Should Know



Tip: Let your system administrator know what you need

Because your PIVOT is deployed within your workplace, it functions differently than a consumer-grade smartphone. It features productivity applications that provide an IT administrator with the ability to restrict any potentially insecure exposure or prevent entertainment-type uses. But these same protections might restrict usage that you need to work smarter.

The Android platform makes it possible to tailor applications that provide functionality specifically designed for your needs.

Therefore, if you find that you need a better, smarter, faster way of conducting your tasks, be sure to alert your system administrator. Chances are that an app is being developed to give you that exact solution.

No camera?

Your PIVOT handset does not have a camera and there is no camera app on the handset. Some apps require some camera-type code and these will not work on PIVOT.

Can't access the Web?

Your system administrator may limit your access to websites. You may be able to browse only to an in house website but not the World Wide Web.

Contact your system administrator if you need information from the World Wide Web.

Can't connect to a computer?

A computer may not be able to recognize your PIVOT handset when you plug your PIVOT into the computer through a USB cable.

Contact your system administrator if you need to be able to access information from a computer on your PIVOT handset.

Can't print?

Although there is a Printing option in the Settings menu, the services that it requires may not be installed. Contact your system administrator.

Appendix A: Care of the Handset

Your PIVOT handset is constructed for durability in a workplace environment. Follow these guidelines to maximize its expected lifetime.



Caution

There are no serviceable parts in the handset or chargers. You should not open the handset case nor disassemble the charger. Doing so will void your warranty.

- Spectralink recommends the use of a silicone bumper to prolong the life of your handset, especially if it could be dropped or cleaned frequently with disinfectants that are known to damage plastics. Try to avoid dropping the handset or knocking it against hard surfaces.
- Turn off the handset and unplug the chargers before you clean them. Never immerse the handset or charger in liquid.
- Take care not to exert undue pressure on the glass screen or on the electrical contacts on the handset, Battery Pack or charger.
- Wiping the handset surface with a water-dampened cloth or paper towel will remove most films or residues. If the soiling is too stubborn for plain water, a mild detergent solution, Lysol, isopropyl alcohol, or diluted bleach may be used on the plastic. Avoid spraying any solution directly onto the handset. Dampen a cloth and wipe instead.
- Pre-treated cloths and towelettes, like those used for eyeglasses or cameras, may be used to clean the glass screen. Avoid using those containing lanolin or aloe as these will leave a slippery residue.
- Use caution when cleaning the Spectralink PIVOT handset with any cleaning products containing quaternary ammonium compounds, such as CaviWipes™, CaviWipes XL, Cavicide® or Sani-Wipes®, Sani-Cloth®, Super Sani-Cloth and Sani-Cloth Plus.
- Do not mix cleaning agents. The combined effects of cleaning agents are unknown. Mixing chemical agents could seriously degrade the construction of the handset making it susceptible to damage, even with normal use.
- Do not use these products to clean your PIVOT handset:
 - Do not use furniture polishes, waxes, or plasticizer-based cleaners (ArmorAll®, etc.)
 - Do not use lanolin, aloe, glycerin, or other skin care products
 - Do not apply any solvent such as acetone, mineral spirits, etc.
- Should the headset connector become dirty, you may experience a scratchy or intermittent signal. To clean the connector, dip the non-padded end of either a wooden or paper handled cotton swab in isopropyl alcohol. Gently insert in the connector and twist, repeating several times. If available, blow compressed air into the connector to clear debris.

- The headset jack and the USB port should be kept covered with the rubber plugs when they are not in use. This will protect the handset against ingress of dust and liquids.
- The handset has an IP64 rating. This means it can withstand dust and splashing liquids but is not submersible. Exposing the unit to liquids may result in permanent damage to the unit.
- If your wireless handset interior gets wet, do not try to accelerate drying with the use of an oven or a dryer as this will damage the wireless handset and void the warranty. Instead, do the following:
 - a** Immediately power off the handset and remove Battery Pack.
 - b** Shake excess liquid from the handset.
 - c** Place the handset and Battery Pack in an area that is at room temperature and has good airflow.
 - d** Let the handset and Battery Pack dry for 72 hours before reconnecting the Battery Pack and/or powering on the handset.

If the handset does not work after following the steps listed above, contact your system administrator.

Appendix B: Regulatory and Safety Information

Product Compatibility



Caution: Product compatibility

Each Spectralink family of products only works with other members of the same family. If you have any questions about product compatibility, contact your system administrator.

Spectralink 87-Series

Use only 87-Series products with other 87-Series products as identified by the model number located on the label of the product.

Spectralink 8753 Wireless Telephone	875X
Spectralink 8742 Wireless Telephone	874X
Spectralink 8741 Wireless Telephone	874X
USB Charger	UE08WCP-050150SPA
Quad Charger	CBK87300 (for both DM351 and DM367)
Quad Charger	CBK87200 (for DM351)
Standard Capacity Battery Pack S:	DM367 (for 8742)
Standard Capacity Battery Pack	DM351 (for 8741/8753)

Spectralink 84-Series

Use only 84-Series products with other 84-Series products as identified by the model number located on the label of the product.

Spectralink 8453 Wireless Handset	8453
Spectralink 8452 Wireless Handset	8452
Spectralink 8450 Wireless Handset	8450
Spectralink 8441 Wireless Handset	8441
Spectralink 8440 Wireless Handset	8440
Speakerphone Dock	8410
USB Charger	SA106B-05
Dual Charger	DCA39
Quad Charger	DCA40
12-Bay Multi-Charger	S0001
Standard Capacity Battery Pack	RS657
Extended Capacity Battery Pack	RS658

Spectralink 80-Series

Use only 80-Series products with other 80-Series products as identified by the model number located on the label of the product.

Spectralink 8020 Wireless Handset	802X
Spectralink 8030 Wireless Handset	803X
Desktop Charger	PCS1850
Dual Charger	PCD1850
Quad Charger	PCQ1850
Standard Battery Pack	PBP0850
Extended Battery Pack	PBP1300
UltraExtended Battery Pack	PBP1850

FCC

FCC IDs IYG-874X and IYG-875X

FCC Radiation Exposure Statement This portable device with its antenna complies with FCC's RF radiation exposure limits set forth for an uncontrolled environment. To maintain compliance this transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

FCC Compliance Statement This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: 1) This device may not cause harmful interference, and 2) this device must accept any interference received, including interference that may cause undesired operation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.



Industry Canada Notice

Certification Number IC: 2128B-874

Certification Number IC: 2128B-875

Under Industry Canada regulations, this radio transmitter may only operate using an antenna of a type and maximum (or lesser) gain approved for the transmitter by Industry Canada. To reduce potential radio interference to other users, the antenna type and its gain should not exceed the equivalent isotropically radiated power (EIRP) necessary for successful communication.

Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes : (1) l'appareil ne doit pas produire de brouillage, et (2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.

This Class B digital apparatus complies with Canadian ICES-003.

Cet appareil numérique de la classe B est conforme à la norme NMB-003 du Canada.

This device and its antenna(s) must not be co-located or operating in conjunction with any other antenna or transmitter.

Cet appareil et son antenne (s) ne doit pas être co-localisés ou fonctionnant en conjonction avec une autre antenne ou transmetteur.

This equipment should be installed and operated with a minimum distance of 20cm between the radiator and your body.

Cet équipement doit être installé et utilisé à une distance minimale de 20cm entre le radiateur et votre corps.

It is the responsibility of the installer to ensure that when using the outdoor antenna kits, only those antennas certified with the product are used. The use of any antenna other than those certified with the product is expressly forbidden by FCC rules 47 CFR part 15.204 and IC RSS standards.

Il est de la responsabilité de l'installateur de s'assurer que lorsque vous utilisez les kits d'antennes extérieures, seules les antennes certifiées avec le produit sont utilisés. L'utilisation d'une antenne autre que ceux qui sont certifiés avec le produit est expressément interdite par la réglementation FCC partie 47 CFR 15.204 et IC normes RSS.

CE Declaration of Conformity

CE Declaration of Conformity with Regard to Directives: 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Spectralink Model No.: 874X, 875X

Български:

С настоящето, Spectralink Corporation декларира, че този 874X, 875X отговаря на основните изисквания и други съответни постановления на Директива 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Čeština:

Spectralink Corporation tímto prohlašuje, že tento 874X, 875X je ve shodě se základními požadavky a dalšími příslušnými ustanoveními směrnice 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Dansk:

Undertegnede Spectralink Corporation erklærer herved, at følgende udstyr 874X, 875X overholder de væsentlige krav og øvrige relevante krav i direktiv 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Deutsch:

Hiermit erklärt Spectralink Corporation, dass sich dieses 874X, 875X in Übereinstimmung mit den grundlegenden Anforderungen und den anderen relevanten Vorschriften der **Richtlinie** 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU befindet.

Eesti:

Käesolevaga kinnitab Spectralink Corporation seadme 874X, 875X vastavust direktiivi 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU põhinõuetele ja nimetatud direktiivist tulenevatele teistele asjakohastele sätetele.

Ελληνικά:

ΜΕ ΤΗΝ ΠΑΡΟΥΣΑ Spectralink Corporation ΔΗΛΩΝΕΙ ΟΤΙ 874X, 875X ΣΥΜΜΟΡΦΩΝΕΤΑΙ ΠΡΟΣ ΤΙΣ ΟΥΣΙΩΔΕΙΣ ΑΠΑΙΤΗΣΕΙΣ ΚΑΙ ΤΙΣ ΛΟΙΠΕΣ ΣΧΕΤΙΚΕΣ ΔΙΑΤΑΞΕΙΣ ΤΗΣ ΟΔΗΓΙΑΣ 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

English:

Hereby, Spectralink Corporation, declares that this 874X, 875X is in compliance with the essential requirements and other relevant provisions of Directive 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Español:

Por medio de la presente Spectralink Corporation declara que el 874X, 875X cumple con los requisitos esenciales y cualesquiera otras disposiciones aplicables o exigibles de la Directiva 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Français:

Par la présente Spectralink Corporation déclare que l'appareil 874X, 875X est conforme aux exigences essentielles et aux autres dispositions pertinentes de la directive 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Gaeilge:

Fograíonn Spectralink Corporation leis seo go bhfuil an 874X, 875X seo i gcomhlíonadh leis na fíor-riachtanais agus na forálacha eile maidir le Treoir 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Hrvatski:

Ovime tvrtka Spectralink Corporation izjavljuje da je ovaj 874X, 875X u skladu osnovnim zahtjevima I ostalim odredbama Direktive 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Íslenska:

Hér með lýsir Spectralink Corporation yfir því að 874X, 875X er í samræmi við grunnkröfur og aðrar kröfur, sem gerðar eru í tilskipun 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Italiano:

Con la presente Spectralink Corporation dichiara che questo 874X, 875X è conforme ai requisiti essenziali ed alle altre disposizioni pertinenti stabilite dalla direttiva 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Latviešu valoda:

Ar šo Spectralink Corporation deklarē, ka 874X, 875X atbilst Direktīvas 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU būtiskajām prasībām un citiem ar to saistītajiem noteikumiem.

Lietuvių kalba:

Šiuo Spectralink Corporation deklaruoja, kad šis 874X, 875X atitinka esminius reikalavimus ir kitas 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU Direktyvos nuostatas.

Magyar:

Alulírott, Spectralink Corporation nyilatkozom, hogy a 874X, 875X megfelel a vonatkozó alapvető követelményeknek és az 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU irányelv egyéb előírásainak.

Malti:

Hawnhekk, Spectralink Corporation jiddikjara li dan 874X, 875X jikkonforma mal-ħtiġijiet essenzjali u ma provvedimenti oħrajn relevanti li hemm fid-Dirrettiva 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Nederlands:

Hierbij verklaart Spectralink Corporation dat het toestel 874X, 875X in overeenstemming is met de essentiële eisen en de andere relevante bepalingen van richtlijn 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Norsk:

Spectralink Corporation erklærer herved at utstyret 874X, 875X er i samsvar med de grunnleggende krav og øvrige relevante krav i direktiv 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Polski:

Niniejszym Spectralink Corporation oświadcza, że 874X, 875X jest zgodny z zasadniczymi wymogami oraz pozostałymi stosownymi postanowieniami Dyrektywy 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Português:

Spectralink Corporation declara que este 874X, 875X está conforme com os requisitos essenciais e outras disposições da Directiva 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Română:

Prin prezenta, Spectralink Corporation declara ca acest 874X, 875X este in conformitate cu cerintele esentiale si alte prevederi ale Directivei 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Slovenčina:

Spectralink Corporation týmto vyhlasuje, že 874X, 875X spĺňa základné požiadavky a všetky príslušné ustanovenia Smernice 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Slovenščina:

Spectralink Corporation izjavlja, da je ta 874X, 875X v skladu z bistvenimi zahtevami in ostalimi relevantnimi določili direktive 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Suomi:

Spectralink Corporation vakuuttaa täten että 874X, 875X tyyppinen laite on direktiivin 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU oleellisten vaatimusten ja sitä koskevien direktiivin muiden ehtojen mukainen.

Svenska:

Härmed intygar Spectralink Corporation att denna 874X, 875X står i överensstämmelse med de väsentliga egenskapskrav och övriga relevanta bestämmelser som framgår av direktiv 2006/95/EC, 2004/108/EC, 1999/5/EC & 2011/65/EU.

Hearing Aid Compatibility

This equipment is Hearing Aid Compatible (HAC)

Cet équipement est compatible avec les dispositifs d'aide auditive.

Chargers and Battery Packs

- Chargers operate in a 50° to 85° F (10° to 30° C) environment. Do not expose them to freezing temperatures or direct sunlight.
- Do not place anything in the charger other than the handset or single Battery Pack as appropriate to avoid bending the contacts. Bent contacts can keep the Battery Pack from charging.
- Only use Spectralink 87-Series Battery Packs with Spectralink 87-Series Chargers.
- Never use a non-Spectralink charging unit as it could damage the Battery Pack.
- Use only the original plug-in power adapter for the chargers.
- Contact your service representative for assistance if you have a faulty Battery Pack or a problem with your charger.
- Any battery that exhibits swelling, cracking or other abnormality should be disposed of promptly and properly. Do not use in such a condition.
- Store Battery in dry conditions at approx. 65° F (20° C).
- Recharge a stored Battery Pack at least every six months.
- Never store a Battery Pack in fully or near fully discharged state.



Warning: Risk of explosion

Use only Spectralink Battery Pack model DM351 with Spectralink 87-Series handsets. Use of another battery may present a risk of fire or explosion. Use Chargers only as specified in this document.



Warning: Do not mistreat

The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not disassemble, crush, heat above 86° F (30° C), or incinerate.



Warning: Water or fire exposure

Do not immerse the Battery Pack in liquid or throw it into fire.



Warning: Disposal of Battery Packs

Do not throw away the Battery Pack with your domestic waste. Return used Battery Packs to your supplier or servicing agent. The battery must be discarded according to Instructions for Battery Collection for each Local Area



Warning: Proper handling

All batteries can cause property damage and/or bodily injury, such as burns, if a conductive material such as jewelry, keys, or beaded chains touches exposed terminals. The conductive material may complete an electrical circuit (short circuit) and become quite hot. Exercise care in handling any charged battery, particularly when placing it inside a pocket, purse, or other container with metal objects.

Wireless Telephones

Please read the following text carefully for important safety information.



Warning: Magnetic earpiece

The earpiece may retain magnetic objects.



Caution: Authority to operate this equipment

Changes or modifications to this equipment that are not approved by Spectralink may cause this equipment to fail to comply with part 15 of the FCC rules, voiding the user's authority to operate this equipment.



Caution: Not user serviceable

Spectralink products contain no user-serviceable parts inside. Refer servicing to qualified service personnel.



Warning: Medical equipment interference

Spectralink recommends that standard acceptance procedures be followed prior to operating this equipment in proximity of life-support systems.

Operational Warnings

Certain operational environments are potentially hazardous. Please read the following text carefully to ensure that you understand the handset's operating conditions.

Potentially Explosive Atmospheres

Do not take your handset into any area with a potentially explosive atmosphere. Do not remove, install, or charge batteries in such areas. Sparks in a potentially explosive atmosphere can cause an explosion or fire resulting in bodily injury or even death.



Warning: Explosive Atmospheres

Avoid areas with potentially explosive atmospheres include fueling areas such as lower decks on boats, fuel, or chemical transfer or storage facilities, areas where the air contains chemicals or particles such as grain, dust, or metal powders, and any other area where you would normally be advised to turn off your vehicle engine. Areas with potentially explosive atmospheres are often marked with signs, but not always.

Electromagnetic Interference and Compatibility

Nearly every electronic device is susceptible to electromagnetic interference (EMI) if inadequately shielded, designed, or otherwise configured for electromagnetic compatibility.

Facilities

To avoid electromagnetic interference and/or compatibility conflicts, turn off your handset in any facility where posted notices instruct you to do so. Hospitals or health care facilities may be using equipment that is sensitive to external radio frequency (RF) energy.

Pacemakers

The Health Industry Manufacturers Association recommends that you maintain a minimum separation of 6 inches (15 cm) between a handheld, wireless-radio product and a pacemaker. These recommendations are consistent with the independent research by, and recommendations of, Wireless Technology Research. Persons with pacemakers should:

- ALWAYS keep the handset more than 6 inches (15 cm) from their pacemaker when the handset is turned on
- Not carry the handset in a breast pocket
- Use the ear opposite the pacemaker to minimize the potential for interference
- Turn the handset off immediately if you have any reason to suspect that interference is taking place

Other Medical Devices

If you use any other personal medical device, consult the manufacturer of your device to determine if it is adequately shielded from external RF energy. Your physician may be able to assist you in obtaining this information.

Use While Driving

Campus installations and warehouse facilities that use vehicles such as forklifts or golf carts should abide by these guidelines when using wireless telephones:

- Give full attention to driving and to the road, aisle, or path
- Use hands-free operation, if available
- Pull off the road, aisle, or path and park before making or answering a call

For Vehicles Equipped with an Airbag

Do not place a portable device in the area over the airbag or in the airbag deployment area. An airbag inflates with great force. If you place the handset in the airbag deployment area, and the airbag inflates, the airbag may propel the handset at high speed and cause serious injury to occupants of the vehicle.

FCC Statement

This device complies with part 15 of the United States FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Specific Absorption Rate (SAR) Information: Spectralink 87-Series Wireless Handsets

Your wireless handset is a low power radio transmitter and receiver. When it is on, it receives and sends out 802.11a/b/g/n radio frequency (RF) signals. In August 1996, the Federal Communications Commission (FCC) adopted radio frequency (RF) exposure guidelines with safety levels for hand-held wireless handsets. Those guidelines are consistent with the safety standards previously set by both U.S. and international standards bodies:

- ANSI C95.1 (1992) American National Standards Institute
- NCRP - Report 86 (1986) National Council on Radiation Protection and Measurements
- ICNIRP (1996) International Commission on Non-Ionizing Radiation Protection
- DHWC - Safety Code 6 Department of Health and Welfare Canada

These standards were developed by independent scientific organizations through periodic and thorough evaluation of scientific studies. The standards include a substantial safety margin

designed to assure the safety of all persons, regardless of age and health. The exposure standard for wireless mobile handsets employs a unit of measurement known as the Specific Absorption Rate, or SAR. The SAR limit set by the FCC is 1.6W/kg.

**Note: SAR limits**

In the United States and Canada, the SAR limit for mobile handsets used by the public is 1.6 watts/kg (W/kg) averaged over one gram of tissue. The standard incorporates a substantial margin of safety to give additional protection for the public and to account for any variations in measurements.

The FCC conducts tests for SAR using standard operating positions specified by the FCC with the handset transmitting at its highest certified power level in all tested frequency bands. Although the SAR is determined at the highest certified power level, the actual SAR level of the handset while operating can be well below the maximum value. This is because the handset is designed to operate at multiple power levels so that it uses only the power required to reach the network. Before a handset model is available for sale to the public, the FCC must test and certify it to ensure that handsets do not exceed the limit established by the government-adopted requirement for safe exposure. The tests are performed in positions and locations (at the ear and worn on the body) as required by the FCC for each model. While there may be differences between the SAR levels of various handsets and at various positions, they all meet the government requirement for safe exposure.

The FCC has granted an Equipment Authorization for this model handset with all reported SAR levels evaluated in compliance with the FCC RF emission guidelines.

SAR information on this model handset is on file with the FCC. You can find it under the Display Grant section at <http://www.fcc.gov/oet/ea> after searching on the FCC ID number found on the handset label inside the battery compartment. Additional information on Specific Absorption Rates (SAR) can be found on the Cellular Telecommunications Industry Association (CTIA) Website at <http://www.ctia.org>.

EU SAR peak values (in W/kg):

8741/53-head: 0.44 W/kg Head Tissue

8741/53-body: 0.48 W/kg Body Tissue

8742-head: 0.51 W/kg Head Tissue

8742-body: 0.60 W/kg Body Tissue

FCC SAR peak values (in W/kg):

8741/53-head: 1.320 W/kg Head Tissue

8741/53-body: 0.280 W/kg Body Tissue

8742-head: 1.420 W/kg Head Tissue

8742-body: 0.410 W/kg Body Tissue

Handset operation normal position

To position the antenna properly, hold the handset as you would any other telephone, with the earpiece to your ear and speak into the microphone.

Appendix C: Products Mentioned in this Document

Android, Google, Google Play and other marks are trademarks of Google, Inc.

Microsoft and Outlook are registered trademarks of Microsoft Corporation in the United States and/or other countries.