

pEF1 α -DsRed-Monomer-N1 Vector

Table of Contents

Product Information	1
Description.....	2
Location of Features	3
Additional Information	3
Quality Control Data.....	4

Catalog No.

631978

Amount10 μ g**Lot Number**

Specified on product label.

Product Information

pEF1 α -DsRed-Monomer-N1 is a mammalian expression vector that constitutively expresses a protein of interest fused to the N-terminus of the red fluorescent protein DsRed-Monomer, even after stable integration of the vector into the host cell genome. Stable, constitutive expression of the fusion protein is driven by the human elongation factor 1 alpha (EF1 α) promoter, allowing the monitoring of a variety of cellular processes (such as differentiation in primary or stem cells) without the transgene silencing associated with CMV promoters. The unmodified vector can be used to express DsRed-Monomer in mammalian cells.

Package Contents

- 1 tube of pEF1 α -DsRed-Monomer-N1 Vector (20 μ l/tube)

Storage Conditions

- Store plasmid at -20°C .
- Spin briefly to recover contents.
- Avoid repeated freeze/thaw cycles.

Shelf Life

- 1 year from date of receipt under proper storage conditions.

Storage Buffer

- 10 mM Tris-HCl (pH 8.0), 1 mM EDTA (pH 8.0)

Concentration

- 500 ng/ μ l

Shipping Conditions

- Dry ice (-70°C)

Clontech Laboratories, Inc.

A Takara Bio Company

1290 Terra Bella Avenue, Mountain View, CA 94043, USA

U.S. Technical Support: tech@clontech.com

United States/Canada

800.662.2566

(PA124336)

Asia Pacific

+1.650.919.7300

Europe

+33.(0)1.3904.6880

Japan

+81.(0)77.543.6116

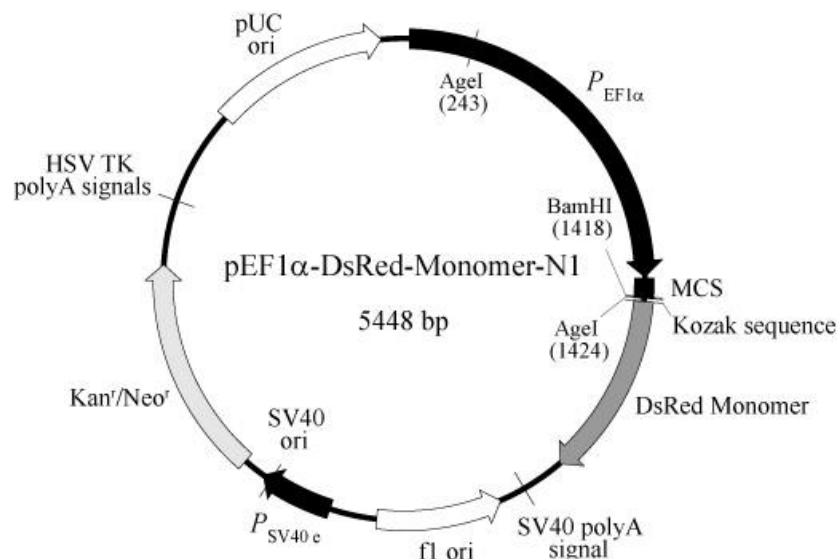


Figure 1. pEF1 α -DsRed-Monomer-N1 vector map. Please note that the vector DNA provided by Clontech is methylated. If you wish to digest the vector with methylation-sensitive enzymes, you will first need to transform the vector into a *dam*⁻ host strain and purify fresh plasmid DNA.

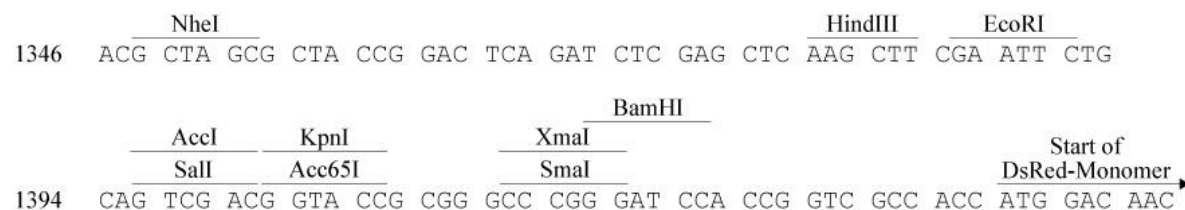


Figure 2. pEF1 α -DsRed-Monomer-N1 multiple cloning site (MCS).

Description

pEF1 α -DsRed-Monomer-N1 is designed to express a protein of interest fused to the N-terminus of DsRed-Monomer, a monomeric mutant of the *Discosoma sp.* red fluorescent protein DsRed (1). The DsRed-Monomer coding sequence has been human-codon-optimized for high expression in mammalian cells (2). The excitation and emission maxima of native DsRed-Monomer are 557 nm and 585 nm, respectively. Expression of fusion proteins that retain the fluorescence properties of unmodified DsRed-Monomer can be monitored by flow cytometry and localized by fluorescence microscopy.

The multiple cloning site (MCS) in pEF1 α -DsRed-Monomer-N1 is positioned between the EF1 α promoter ($P_{EF1\alpha}$) and the DsRed-Monomer coding sequence. Expression of the fusion protein is driven by the EF1 α promoter, which remains constitutively active even after stable integration of the vector into the host cell genome (3). A Kozak consensus sequence, located immediately upstream of the DsRed-Monomer gene, enhances the translational efficiency of DsRed-Monomer in eukaryotic systems (4), and SV40 polyadenylation signals direct proper processing of the 3' end of the DsRed-Monomer mRNA.

The vector backbone contains an SV40 origin for replication in mammalian cells expressing the SV40 large T antigen, a pUC origin of replication for propagation in *E. coli*, and an f1 origin for single-stranded DNA production. A neomycin-resistance cassette (Neo^r) allows stably transfected eukaryotic cells to be selected using G418 (5). This cassette consists of

the SV40 early promoter ($P_{SV40\ e}$), the Tn5 neomycin/kanamycin resistance gene, and polyadenylation signals from the herpes simplex virus thymidine kinase (HSV TK) gene. A bacterial promoter upstream of the cassette drives expression of the kanamycin resistance gene in *E. coli*.

Location of Features

- $P_{EF1\alpha}$ (human elongation factor 1 alpha promoter): 12–1346
- MCS (multiple cloning site): 1348–1422
- Kozak consensus sequence: 1429–1439
- DsRed-Monomer (human-codon-optimized): 1436–2110
- SV40 polyA signal: 2267–2301
- f1 origin of replication: 2364–2819 (complementary)
- $P_{SV40\ e}$ (SV40 early promoter and enhancer sequences): 2993–3261
- SV40 origin of replication: 3160–3295
- Kan^r/Neo^r (kanamycin/neomycin resistance gene): 3344–4138
- HSV TK polyA signals: 4374–4392
- pUC origin of replication: 4723–5366

Additional Information

The gene of interest must be cloned into pEF1 α -DsRed-Monomer-N1 so that it is in-frame with the DsRed-Monomer coding sequence. The gene must contain a start codon (ATG), and lack in-frame stop codons.

Cells expressing DsRed-Monomer fusions can be detected by flow cytometry or fluorescence microscopy 12–16 hr after transfection. If required, stable transfectants can be selected using G418 (5). pEF1 α -DsRed-Monomer-N1 can also be used as a cotransfection marker, as the unmodified vector will express DsRed-Monomer in mammalian cells.

Propagation in *E. coli*

- Suitable host strains: DH5 α , HB101 and other general purpose strains. Single-stranded DNA production requires a host containing an F plasmid, such as the JM109 or XL1-Blue strains.
- Selectable marker: plasmid confers resistance to kanamycin (50 μ g/ml) in *E. coli* hosts.
- *E. coli* replication origin: pUC
- Copy number: high

Excitation and Emission Maxima of DsRed-Monomer

- Excitation: 557 nm
- Emission: 585 nm

References

1. Matz, M. V., *et al.* (1999) *Nat. Biotechnol.* **17**(10):969–973.
2. Haas, J., *et al.* (1996) *Curr. Biol.* **6**(3):315–324.
3. Wang, R. *et al.* (2008) *Stem Cells Dev.* **17**(2):279–289.
4. Kozak, M. (1987) *Nucleic Acids Res.* **15**(20):8125–8148.
5. Gorman, C. (1985) In *DNA Cloning: A Practical Approach, Vol. II*. Ed. D. M. Glover (IRL Press, Oxford, U.K.) pp. 143–190.

Certificate of Analysis

Cat. No. 631978

pEF1 α -DsRed-Monomer-N1 Vector

Quality Control Data

Plasmid Identity & Purity

- Digestion with the indicated restriction enzymes produced fragments of the indicated sizes on a 0.8% agarose/EtBr gel:

Enzyme(s)	Fragment(s)
BamHI	5.4 kb
AgeI	1.2 & 4.3 kb

- Vector identity was confirmed by sequencing.
- A_{260}/A_{280} : 1.8–2.0

pEF1alpha-DsRed-Monomer-N1 Vector

CATALOG NO.

631978

NOTICE TO PURCHASER:

Clontech products are to be used for research purposes only. They may not be used for any other purpose, including, but not limited to, use in drugs, in vitro diagnostic purposes, therapeutics, or in humans. Clontech products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products or to provide a service to third parties without prior written approval of Clontech Laboratories, Inc.

Your use of this product is also subject to compliance with the licensing requirements listed below and described on the product's web page at <http://www.clontech.com>. It is your responsibility to review, understand and adhere to any restrictions imposed by these statements.

STATEMENT 44

The DsRed-Monomer and the Fruit Fluorescent Proteins are covered by one or more of the following U.S. Patents: 7,005,511; 7,157,566; 7,393,923 and 7,250,298.

STATEMENT 72

Living Colors Fluorescent Protein Products: Not-For-Profit Entities: Orders may be placed in the normal manner by contacting your local representative or Clontech Customer Service at. Any and all uses of this product will be subject to the terms and conditions of the Non-Commercial Use License Agreement (the "Non-Commercial License"), a copy of which can be found below. As a condition of sale of this product to you, and prior to using this product, you must agree to the terms and conditions of the Non-Commercial License. Under the Non-Commercial License, Clontech grants Not-For-Profit Entities a non-exclusive, non-transferable, non-sublicensable and limited license to use this product for internal, non-commercial scientific research use only. Such license specifically excludes the right to sell or otherwise transfer this product, its components or derivatives thereof to third parties. No modifications to the product may be made without express written permission from Clontech. Any other use of this product requires a different license from Clontech. For license information, please contact a licensing representative by phone at 650.919.7320 or by e-mail at licensing@clontech.com. For-Profit Entities wishing to use this product are required to obtain a license from Clontech. For license information, please contact a licensing representative by phone at 650.919.7320 or by e-mail at licensing@clontech.com.

TRADEMARKS:

Clontech Laboratories, Inc.

A Takara Bio Company

1290 Terra Bella Avenue, Mountain View, CA 94043, USA

U.S. Technical Support: tech@clontech.com

4/3/2014

United States/Canada

800.662.2566

Asia Pacific

+1.650.919.7300

Europe

+33.(0)1.3904.6880

Japan

+81.(0)77.543.6116

Notice to Purchaser



Clontech and the Clontech logo are trademarks of Clontech Laboratories, Inc.

All other marks are the property of their respective owners. Certain trademarks may not be registered in all jurisdictions. Clontech is a Takara Bio Company. ©2014 Clontech Laboratories, Inc. This document has been reviewed and approved by the Clontech Quality Assurance Department.

Clontech Laboratories, Inc.

A Takara Bio Company

1290 Terra Bella Avenue, Mountain View, CA 94043, USA

U.S. Technical Support: tech@clontech.com

4/3/2014

United States/Canada

800.662.2566

Asia Pacific

+1.650.919.7300

Europe

+33.(0)1.3904.6880

Japan

+81.(0)77.543.6116