



5

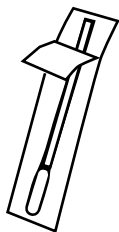

Remove the swab from the package and do not touch the tip of the swab! Open the specimen tube, holding it in one hand and the swab in the other hand. Do not spill or drink the liquid.

6

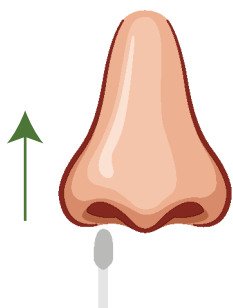


Insert the swab into one nostril.

- For adults - insert less than **one** inch.
- For younger children - insert less than 1/2 inch.

Rotate the swab in a circular motion around the entire inside edge of the nostril. Do this 2 times and then keep it in place for 15 seconds.

Repeat the same process in the other nostril using the same swab.

7

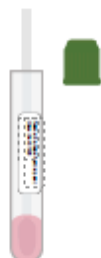


Insert the swab into the tube (tip first) until it reaches the bottom. Tip of swab will be in the liquid. Break the swab handle against the side of the tube at perforated break point.

8



Screw the cap onto the tube and place the tube into the bag. Place the completed requisition into the back pocket of the bag, separate from the tube. The test requisition **MUST** be returned with the specimen in order to be tested. Seal the bag. **Important:** Do not remove the desiccant (white sheet) from the bag.

9

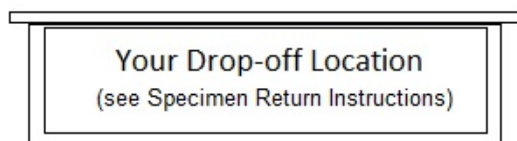

Drop-off the specimen at a drop-off location ideally the same day that the specimen was collected - but no later than 24 hours after collecting the specimen. DO NOT drop-off the specimen if the drop-off location is closed. See the Specimen Return Instructions provided by the facilitator of the testing program.



Collection Kit for COVID-19 (Bulk)

Read all instructions before starting specimen collection.

Rx Only, IVD, and For Use Under Emergency Use Authorization

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of medical devices under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

LabOne, Inc. dba Quest Diagnostics, Lenexa KS, Quest Diagnostics Incorporated and its subsidiaries (Quest) complies with applicable federal civil rights laws and does not discriminate on the basis of race, color, national origin, age, disability, or sex. ATTENTION: If you speak English, language assistance services, free of charge, are available to you. ATENCIÓN: Si habla español (Spanish), tiene a su disposición servicios gratuitos de asistencia lingüística. 注意：如果您使用繁體中文 (Chinese), 您可以免費獲得語言援助服務。QuestDiagnostics.com/home/nondiscrimination/

Quest, Quest Diagnostics, any associated logos, and all associated Quest Diagnostics registered or unregistered trademarks are the property of Quest Diagnostics. All third-party marks-® and ™- are the property of their respective owners. ©2021 Quest Diagnostics Incorporated. All rights reserved.

INS0002 10/8/2021

Welcome to Specimen Collection

Before starting to collect a specimen, view the specimen collection demo video at

Patient.QuestDiagnostics.com/COVID-2



READ FIRST FOR YOUR SAFETY

- For use by individuals age 18 years or older (self-collected), 16 years and older (self-collected under adult supervision, i.e., individuals 16-17 years old and adults that require supervision), or 2 years and older (collected with adult assistance, i.e., individuals 2-15 years old and adults that require assistance).
- Do not use if the individual being swabbed is on blood thinners or had a recent nose injury.
- Read the instructions carefully before starting.
- Once the swab is inserted, the nose may feel very full. If there is pain at any time during the collection, stop and call the doctor to guide the specimen collection.
- Do not use the liquid in the tube to wet the swab, or to lubricate the nostrils. Do not drink the liquid.
- Do not touch the swab tip with anything except the inside of the nose and the tube liquid. If the swab tip touches anything, then contact Quest Diagnostics at 1.855.332.2533 for another swab.
- Do not use the swab on more than one person. Sharing a swab with someone else can cause an infection.

Read all instructions before starting specimen collection.

Failure to follow instructions may lead to incorrect test results.

Before you begin

See the Specimen Return Instructions provided by the facilitator of your testing to view drop-off locations. Be sure to bring the specimen to a drop-off location on the same day it is collected, and before the location closes. **DO NOT** drop-off the sample if the drop-off location is closed. If you have any questions, please call 1.855.332.2533.

Materials



Test requisition (pre-printed)



Bag (containing a desiccant)



Tube (containing liquid)



Swab (in a wrapper)



Specimen Return Instructions (provided by the facilitator of your testing program)

This product has:

- not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA.
- has been authorized only for the collection and maintenance of anterior nasal swab specimens as an aid in detection of nucleic acid from SARS-CoV2, not for any other viruses or pathogens.

Read all instructions before starting.

- 1 Open the materials, and lay all the materials on a clean surface.

- 2 Locate and sign the test requisition. Verify the patient's name and birth date are correct, and call 1.855.332.2533 if anything is incorrect.

IMPORTANT: This is required information. Complete the Date Collected and Time Collected lines on the test requisition. The specimen may be rejected if the date and time on the test requisition is not completed.




- 3 Wash and dry your hands.

- 4 Peel off the label with the patient's name and date of birth from the bottom of the test requisition and place it on the tube.


IMPORTANT: This is required information.

4

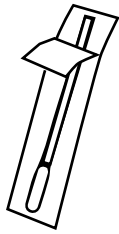

Remove the swab from the package and do not touch the tip of the swab! Open the specimen tube, holding it in one hand and the swab in the other hand. Do not spill or drink the liquid.

5

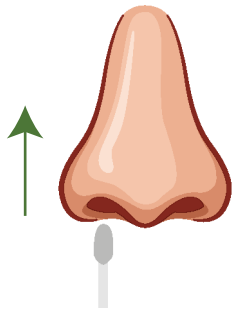

Insert the swab into one nostril.

- For adults - insert less than **one** inch.
- For younger children - insert less than 1/2 inch.

Rotate the swab in a circular motion around the entire inside edge of the nostril. Do this 2 times and then keep it in place for 15 seconds.

Repeat the same process in the other nostril using the same swab.

6



Insert the swab into the tube (tip first) until it reaches the bottom. Tip of swab will be in the liquid. Break the swab handle against the side of the tube at perforated break point.

7



Screw the cap onto the tube and place the tube into the bag. **Important:** Do not remove the desiccant (white sheet) from the bag. Seal the bag.

8



Drop-off the specimen at a drop-off location ideally the same day that the specimen was collected - but no later than 24 hours after collecting the specimen. DO NOT drop-off the specimen if the drop-off location is closed. See the Specimen Return Instructions provided by the facilitator of the testing program.



Collection Kit for COVID-19 (Bulk & Activate)

Read all instructions before starting specimen collection.

Rx Only, IVD, and For Use Under Emergency Use Authorization

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of medical devices under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

LabOne, Inc. dba Quest Diagnostics, Lenexa KS, Quest Diagnostics Incorporated and its subsidiaries (Quest) complies with applicable federal civil rights laws and does not discriminate on the basis of race, color, national origin, age, disability, or sex. ATTENTION: If you speak English, language assistance services, free of charge, are available to you. ATENCIÓN: Si habla español (Spanish), tiene a su disposición servicios gratuitos de asistencia lingüística. 注意: 如果您使用繁體中文 (Chinese), 您可以免費獲得語言援助服務. [QuestDiagnostics.com/home/nondiscrimination/](https://www.questdiagnostics.com/home/nondiscrimination/)

Quest, Quest Diagnostics, any associated logos, and all associated Quest Diagnostics registered or unregistered trademarks are the property of Quest Diagnostics. All third-party marks-® and ™- are the property of their respective owners. ©2021 Quest Diagnostics Incorporated. All rights reserved.

INS0004 10/8/2021

Welcome to Specimen Collection

Before starting to collect a specimen, view the specimen collection demo video at



Patient.QuestDiagnostics.com/COVID-4

READ FIRST FOR YOUR SAFETY

- For use by individuals age 18 years or older (self-collected), 16 years and older (self-collected under adult supervision, i.e., individuals 16-17 years old and adults that require supervision), or 2 years and older (collected with adult assistance, i.e., individuals 2-15 years old and adults that require assistance).
- Do not use if the individual being swabbed is on blood thinners or had a recent nose injury.
- Read the instructions carefully before starting.
- Once the swab is inserted, the nose may feel very full. If there is pain at any time during the collection, stop and call the doctor to guide the specimen collection.
- Do not use the liquid in the tube to wet the swab, or to lubricate the nostrils. Do not drink the liquid.
- Do not touch the swab tip with anything except the inside of the nose and the tube liquid. If the swab tip touches anything, then contact Quest Diagnostics at 1.855.332.2533 for another swab.
- Do not use the swab on more than one person. Sharing a swab with someone else can cause an infection.

Read all instructions before starting specimen collection.

Failure to follow instructions may lead to incorrect test results.

Before you begin

See the Specimen Return Instructions provided by the facilitator of the testing program to view drop-off locations. Be sure to bring the specimen to a drop-off location on the same day it is collected, and before the location closes. **DO NOT** drop-off the sample if the drop-off location is closed. **If you have any questions, please call 1.855.332.2533.**

Materials



Activation Card



Swab (in a wrapper)



Bag (containing a desiccant)



Specimen Return Instructions (provided by the facilitator of your testing program)




Tube (containing liquid)

This product has:

- not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA.
- has been authorized only for the collection and maintenance of anterior nasal swab specimens as an aid in detection of nucleic acid from SARS-CoV2, not for any other viruses or pathogens.

Read all instructions before starting.

- 1 Open the materials, and lay all the materials on a clean surface.

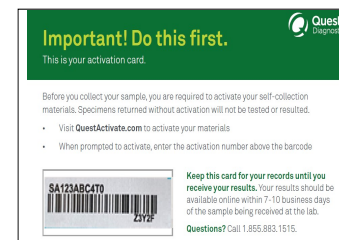

- 2 Locate the Activation Card and follow the instructions on the Activation card to activate your materials the day the specimen is collected. Make sure the code on the Activation Card and the Tube are the same, and call 1.855.332.2533 if the codes are not the same.


IMPORTANT: You must activate the Materials online or the lab cannot test the specimen.

- 3




Wash and dry your hands.

5

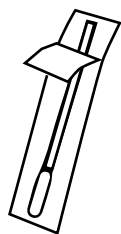

Remove the swab from the package and do not touch the tip of the swab! Open the specimen tube, holding it in one hand and the swab in the other hand. Do not spill or drink the liquid.

6

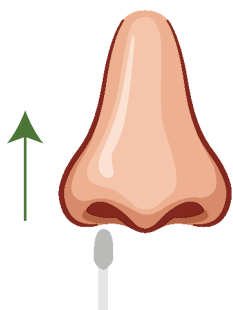


Insert the swab into one nostril.

- For adults - insert less than **one** inch.
 - For younger children - insert less than 1/2 inch.
- Rotate the swab in a circular motion around the entire inside edge of the nostril. Do this 2 times and then keep it in place for 15 seconds.

Repeat the same process in the other nostril using the same swab.

7

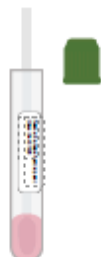

Insert the swab into the tube (tip first) until it reaches the bottom. Tip of swab will be in the liquid. Break the swab handle against the side of the tube at perforated break point.

8



Screw the cap onto the tube and place the tube into the bag. Place the completed requisition into the back pocket of the bag, separate from the tube. The test requisition **MUST** be returned with the specimen in order to be tested. Seal the bag. **Important:** Do not remove the desiccant (white sheet) from the bag. Seal the bag.

9



Place the bag into the box and close the box. Place the box into the Return UN3373 Shipping Bag with the pre-applied FedEx return label. Seal the Shipping Bag. If you are under 18 years of age, have an adult help with packaging and sending the specimen.



Drop the postage-paid FedEx Return Shipping Bag at a drop box location on the same day as the specimen is collected, and before 3:00 PM. **DO NOT** drop the specimen in a drop box on a Saturday or Sunday.



Collection Kit for COVID-19 (Individual)

Read all instructions before starting specimen collection.

Rx Only, IVD, and For Use Under Emergency Use Authorization

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of medical devices under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

LabOne, Inc. dba Quest Diagnostics, Lenexa KS, Quest Diagnostics Incorporated and its subsidiaries (Quest) complies with applicable federal civil rights laws and does not discriminate on the basis of race, color, national origin, age, disability, or sex. ATTENTION: If you speak English, language assistance services, free of charge, are available to you. ATENCIÓN: Si habla español (Spanish), tiene a su disposición servicios gratuitos de asistencia lingüística. 注意：如果您使用繁體中文 (Chinese), 您可以免費獲得語言援助服務。QuestDiagnostics.com/home/nondiscrimination/

Quest, Quest Diagnostics, any associated logos, and all associated Quest Diagnostics registered or unregistered trademarks are the property of Quest Diagnostics. All third-party marks-® and ™- are the property of their respective owners. ©2021 Quest Diagnostics Incorporated. All rights reserved.

INS0001 10/8/2021

Welcome to Specimen Collection

Before starting to collect a specimen, view the specimen collection demo video at



Patient.QuestDiagnostics.com/COVID-1

READ FIRST FOR YOUR SAFETY

- For use by individuals age 18 years or older (self-collected), 16 years and older (self-collected under adult supervision, i.e., individuals 16-17 years old and adults that require supervision), or 2 years and older (collected with adult assistance, i.e., individuals 2-15 years old and adults that require assistance).
- Do not use if the individual being swabbed is on blood thinners or had a recent nose injury.
- Read the instructions carefully before starting.
- Once the swab is inserted, the nose may feel very full. If there is pain at any time during the collection, stop and call the doctor to guide the specimen collection.
- Do not use the liquid in the tube to wet the swab, or to lubricate the nostrils. Do not drink the liquid. Do not touch the swab tip with anything except the inside of the nose and the tube liquid. If the swab tip touches anything, then contact Quest Diagnostics at 1.855.332.2533 for another swab.
- Do not use the swab on more than one person. Sharing a swab with someone else can cause an infection.

Read all instructions before starting specimen collection.

Failure to follow instructions may lead to incorrect test results.

Before you begin

Visit [fedex.com/labreturns](https://www.fedex.com/labreturns) to view drop box locations. Be sure to bring the specimen to a drop box on the same day that the specimen is collected and before 3:00PM. DO NOT drop specimen into a drop box on Saturday or Sunday. If you have any questions, please call 1.855.332.2533.

Materials



Test requisition (pre-printed)



Swab (in a wrapper)



Bag (containing a desiccant)



Tube (containing liquid)



Box



Bag (front)



Bag (back)

Return UN3373 Shipping Bag with pre-applied FedEx shipping label

This product has:

- not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA.
- has been authorized only for the collection and maintenance of anterior nasal swab specimens as an aid in detection of nucleic acid from SARS-CoV2, not for any other viruses or pathogens.

Read all instructions before starting.


- 1 Open the materials, and lay all the materials on a clean surface.

- 2 Locate and sign the test requisition. Verify the patient's name and birth date are correct, and call 1.855.332.2533 if anything is incorrect.

IMPORTANT: This is required information. Complete the Date Collected and Time Collected lines on the test requisition. The specimen may be rejected if the date and time on the test requisition is not completed.




- 3 Wash and dry your hands.

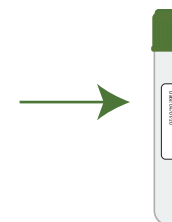

- 4 Peel off the label with the patient's name and date of birth from the bottom of the test requisition and place it on the tube.


IMPORTANT: This is required information.

4

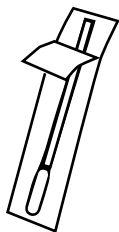

Remove the swab from the package and do not touch the tip of the swab! Open the specimen tube, holding it in one hand and the swab in the other hand. Do not spill or drink the liquid.

5

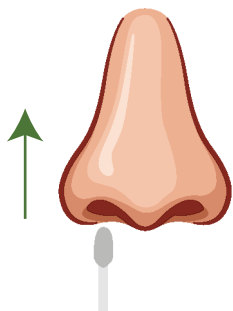

Insert the swab into one nostril.

- For adults - insert less than **one** inch.
- For younger children - insert less than 1/2 inch.

Rotate the swab in a circular motion around the entire inside edge of the nostril. Do this 2 times and then keep it in place for 15 seconds.

Repeat the same process in the other nostril using the same swab.

6



Insert the swab into the tube (tip first) until it reaches the bottom. Tip of swab will be in the liquid. Break the swab handle against the side of the tube at perforated break point.

7



Screw the cap onto the tube and place the tube into the bag. **Important:** Do not remove the desiccant (white sheet) from the bag. Seal the bag.

8



Place the bag into the box and close the box. Place the box into the Return UN3373 Shipping Bag with the pre-applied FedEx return label. Seal the Shipping Bag. If you are under 18 years of age, have an adult help you with packaging and sending the specimen.

9



Drop-off the postage-paid FedEx Return Shipping Bag at a drop box location the same the specimen is collected, and before 3:00PM. DO NOT drop-off the specimen in a drop box on a Saturday or Sunday.



Collection Kit for COVID-19 (Individual & Activate)

**Read all instructions
before starting
specimen collection.**

Rx Only, IVD, and For Use Under Emergency Use Authorization

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of medical devices under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

LabOne, Inc. dba Quest Diagnostics, Lenexa KS, Quest Diagnostics Incorporated and its subsidiaries (Quest) complies with applicable federal civil rights laws and does not discriminate on the basis of race, color, national origin, age, disability, or sex. ATTENTION: If you speak English, language assistance services, free of charge, are available to you. ATENCIÓN: Si habla español (Spanish), tiene a su disposición servicios gratuitos de asistencia lingüística. 注意: 如果您使用繁體中文 (Chinese), 您可以免費獲得語言援助服務. [QuestDiagnostics.com/home/nondiscrimination/](https://www.questdiagnostics.com/home/nondiscrimination/)

Quest, Quest Diagnostics, any associated logos, and all associated Quest Diagnostics registered or unregistered trademarks are the property of Quest Diagnostics. All third-party marks-® and ™- are the property of their respective owners. ©2021 Quest Diagnostics Incorporated. All rights reserved.

INS0003 10/8/2021

Welcome to Specimen Collection

Before starting to collect a specimen, view the specimen collection demo video at



Patient.QuestDiagnostics.com/COVID-3

READ FIRST FOR YOUR SAFETY

- For use by individuals age 18 years or older (self-collected), 16 years and older (self-collected under adult supervision, i.e., individuals 16-17 years old and adults that require supervision), or 2 years and older (collected with adult assistance, i.e., individuals 2-15 years old and adults that require assistance).
- Do not use if the individual being swabbed is on blood thinners or had a recent nose injury.
- Read the instructions carefully before starting.
- Once the swab is inserted, the nose may feel very full. If there is pain at any time during the collection, stop and call the doctor to guide the specimen collection.
- Do not use the liquid in the tube to wet the swab, or to lubricate the nostrils. Do not drink the liquid.
- Do not touch the swab tip with anything except the inside of the nose and the tube liquid. If the swab tip touches anything, then contact Quest Diagnostics at 1.855.332.2533 for another swab.
- Do not use the swab on more than one person. Sharing a swab with someone else can cause an infection.

Read all instructions before starting specimen collection.

Failure to follow instructions may lead to incorrect test results.

Before you begin

Visit fedex.com/labreturns to view drop box locations. Be sure to bring the specimen to a drop box on the same day the specimen is collected and before 3:00 PM. DO NOT drop the specimen into a drop box on Saturday or Sunday.

If you have any questions, please call 1.855.332.2533.

Materials



Activation Card



Bag (containing a desiccant)



Tube (containing liquid)



Swab (in a wrapper)



Box



Bag (front)



Bag (back)


Return UN3373 Shipping Bag with pre-applied FedEx shipping label

This product has:

- not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA.
- has been authorized only for the collection and maintenance of anterior nasal swab specimens as an aid in detection of nucleic acid from SARS-CoV2, not for any other viruses or pathogens.

Read all instructions before starting.

- 1 Open the materials, and lay all the materials on a clean surface.

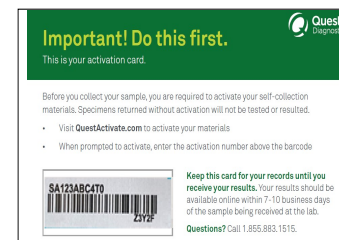

- 2 Locate the Activation Card and follow the instructions on the Activation card to activate your materials the day the specimen is collected. Make sure the code on the Activation Card and the Tube are the same, and call 1.855.332.2533 if the codes are not the same.

IMPORTANT: You must activate the Materials on-line or the lab cannot test the specimen.



- 3 Wash and dry your hands.

Important! Do this first.

This is your activation card.

Before you collect your sample, you are required to activate your collection materials. Specimens returned without activation will not be tested or resulted.

- Visit **QuestActivate.com** to activate your materials
- When prompted to activate, enter the activation number above the barcode



Keep this card for your records until you receive your results. Your results should be available online within 7-10 business days of the sample being received at the lab.

Questions? Call 1.855.883.1515.