

**QUICKIE**<sup>®</sup>

Power Portfolio





## Table of Contents

- Zippie® ZM-310 .....2
- Quickie® QM-7 Series .....3
- Quickie® Pulse Series .....6
- Quickie® S-6 Series .....8
- Quickie® P-222 SE .....10
- PG Drives Technology .....12
- Back Upholstery Options .....12
- Electronics .....13
- Base Options .....16
- Accessories .....17
- Armrests .....18
- Positioning Options .....18
- Legrests .....19
- Built-4-Me .....20
- JAY® Seating Selection .....20

### Why Quickie?

**You drive our design.** Quickie wheelchairs are designed for you and by you. In development, we consider how you interact with your wheelchair to maintain maximum efficiency, performance, fit, and style. This results in the creation of products that are innovative, durable, versatile, and customizable. That is the Quickie Advantage.

**Innovation.** We are continually developing new products, improving our flagship products, and offering new options and solutions to enhance your lifestyle.

**Quality and durability.** Quickie wheelchairs undergo rigorous testing that surpasses industry standards, so you can be confident that your product will hold up and perform throughout its lifetime.

**Versatility.** We offer a comprehensive portfolio of products to meet a wide range of individual lifestyles and clinical needs. Once you choose a Quickie, our goal is that it continues to be the right chair for you. We offer maximum interchangeable parts and adjustability in our designs to help accommodate changes in condition, environment, or lifestyle that you may experience.

**Built-4-Me.** Our Built-4-Me Team is dedicated to creating product customizations beyond what is listed on our product order forms. Experts in the engineering and manufacturing of custom rehab products will listen to your specific requests and build a chair as individual as your fingerprint.

## zippie ZM-310™

The Zippie ZM-310 combines our smallest, most reliable power wheelchair base with clinically proven pediatric seating and positioning options to provide children maximum accessibility and independence.



### Compact Base with Serious Power!

- 24" wide, 35" long base for indoor accessibility
- 19" turning radius for intuitive driving and maneuverability
- 13" drive wheels for climbing obstacles and gliding over terrain

### Positioning Options

- 50° of power tilt
- 50° of manual tilt
- 40° of manual recline
- 9" or 12" of power seat elevate



### Pediatric Seating

- Available in sizes ranging from 10" x 10" to 16" x 19"
- ASAP II seat frame offers in-field width, depth, armrest, legrest, and back height adjustability
- Angle-adjustable centermount footrest with pediatric footplate and 3" to 13" lower leg lengths

### ZIPPIE ZM-310 SPECIFICATIONS

<b>HCPCS Code</b>	K0848, K0856
<b>Drive Wheel</b>	Mid Wheel
<b>Max Speed</b>	6 mph
<b>User Weight Capacity</b>	165 lbs
<b>Width</b>	10" to 16"
<b>Depth</b>	10" to 19"
<b>Seat Height</b>	15.5" to 21.8"
<b>Overall Base Width</b>	24"
<b>Overall Length</b>	35"
<b>Turning Radius</b>	19"
<b>Drive Wheels</b>	13" with flat-free insert
<b>Transit Tested</b>	WC-19

## Live Without Limits

The Quickie QM-7 series offers a high-performing, reliable, and stylish mid-wheel drive power solution for a client who qualifies for a Group 3 Standard (QM-710), Group 3 Heavy Duty (QM-715 HD) or Group 4 High Speed (QM-720) power mobility base.

# QUICKIE® QM-7 Series



## Available in Nine Group 3 and Group 4 Models:

### GROUP 3: Standard Wt. Cap.

- Quickie QM-710 BC (K0848):  
No Power Option Base with ASAP II Rehab Seat
- Quickie QM-710 SC (K0856):  
Single Power Option Tilt (SPOT) or Recline Base with ASAP II or Recline Rehab Seat
- Quickie QM-710 MPC (K0861):  
Multiple Power Option Base (SPOT Tilt and Recline) with 3.7 Recline Seat

### GROUP 3: Heavy Duty

- Quickie QM-715 HD BC (K0850):  
No Power Option Base with ASAP II Rehab Seat
- Quickie QM-715 HD SC (K0858):  
Single Power Option Tilt (SPOT) with ASAP II Rehab Seat
- Quickie QM-715 HD MPC (K0862):  
Multiple Power Option Base (SPOT Tilt and Recline) with 3.7 Recline Seat

### GROUP 4: High Activity

- Quickie QM-720 BD (K0868):  
No Power Option Base with ASAP II Rehab Seat
- Quickie QM-720 SD (K0877):  
Single Power Option Tilt (SPOT) or Recline Base with ASAP II or Recline Rehab Seat
- Quickie QM-720 MPD (K0884):  
Multiple Power Option Base (SPOT Tilt and Recline) with 3.7 Recline Seat

# QUICKIE® QM-7 Series



## Power Base Performance

- Linix 5 mph, 6 mph, or 8 mph, 4-pole motor packages
- Intuitive mid-wheel drive handling
- Contoured, sleek base available in 16 shroud color options
- SpiderTrac™ Suspension for superior traction and handling across multiple terrains
  - Gas dampening system for smooth transitions
  - 3" curb climbing capability
  - 22° transition angle helps prevent high centering on curb cuts



## Seating Options

### ASAP II Seat Frame

- In-field width, depth, armrest, legrest, and back height adjustability
- Wide range of seat sizes
- Manual reclining backrest for increased function and sitting tolerance
- Optional fold down backrest for easy transport
- Extensive seating and positioning options

### 3.7 Recline Seat Frame

- 14" – 22" seat widths and depths
- Ultra low seat-to-floor heights – as low as 16" when combined with tilt and elevate!
- Integrated JAY J3 back for seamless functionality with power recline and shear reduction system
- Cantilever Reclining Armrests for optimal positioning at any back angle
- Low profile electronics enclosure shell for a clean look and compatibility with accessories



## Power Options

- NEW 3.7 Power Recline with Power Shear Reduction
- 50° of Power Tilt to optimize positioning
- 12" and 9" Power Seat Elevate for increased access to cabinets and interaction with peers
- Power Elevating Legrests with Articulation
  - Swing-Away
  - Centermount with Single/Dual Footplate for angle adjustability as a single footplate or separate footplates



## NEW! Options for Consumer Convenience

- **RAM® Phone and Tablet Mounts** allow mobile phones and/or tablets to be connected and disconnected for easy viewing
- **Active Controls Power Buddy** provides an auxiliary USB power source for charging devices and optional LED lights for improved visibility
- **RAM® Mounts Self-Leveling Cup Holder** uses gimbal system to keep drinks level during travel or when stopped on inclines
- **Active Controls ReCharge™ with Socket** eliminates the challenges of connecting and disconnecting the battery charger by providing an easy-to-reach adapter

## Quickie QM-7 Series Specifications

SPECIFICATIONS	QM-710 BC	QM-710 SC	QM-710 MPC	QM-715 HD BC	QM-715 HD SC	QM-715 HD MPC	QM-720 BD	QM-720 SD	QM-720 MPD
HCPCS Code	K0848	K0856	K0861	K0850	K0858	K0862	K0868	K0877	K0884
Drive Wheel Position	Mid	Mid	Mid	Mid	Mid	Mid	Mid	Mid	Mid
Speed	up to 6.0 mph*	up to 6.0 mph*	up to 6.0 mph*	5.0 mph	5.0 mph	5.0 mph	8.0 mph	8.0 mph	8.0 mph
Overall Weight	310 lbs.	350 lbs.	375 lbs.	310 lbs.	350 lbs.	375 lbs.	310 lbs.	350 lbs.	375 lbs.
Overall Width	25"	25"	25"	25"	25"	25"	25"	25"	25"
Overall Length	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests
User Weight Capacity	300 lbs.	300 lbs.	300 lbs.	400 lbs.	Tilt: 400 lbs. Recline: 350 lbs	350 lbs.	300 lbs.	300 lbs.	300 lbs.
Seat Widths	12" to 22"	12" to 22"	14" to 22"	16" to 24"	16" to 24"	16" to 22"	12" to 22"	12" to 22"	14" to 22"
Seat Depths	12" to 22"	12" to 22"	16" to 22"	16" to 22"	16" to 22"	16" to 22"	12" to 22"	12" to 22"	16" to 22"
Seat To Floor Height	16.3" to 21"	16.1" to 20.6"	16" to 19.3"	16.3" to 20.3"	16.1" to 20.8"	16.6" to 19.6"	16.3" to 21"	16.2" to 20.6"	16" to 19.3"
Seat Frame Options	ASAP II	ASAP II or 3.7 Recline	3.7 Recline	ASAP II	ASAP II or 3.7 Recline	3.7 Recline	ASAP II	ASAP II or 3.7 Recline	3.7 Recline
Power Seating Options	N/A	Tilt or Recline	Tilt & Recline	No Power Options	Tilt or Recline	Tilt & Recline	No Power Options	Tilt or Recline	Tilt & Recline

\* Optional 5.0 mph and 8.0 mph motor packages also available

## *Live Without Limits*

The Quickie Pulse is a compact, durable power wheelchair that offers a wide range of effective seating and electronics options to meet the needs of Group 3 users.

# QUICKIE® Pulse™



### **Available in Six Models:**

- Quickie Pulse 5 BC (K0848): 5 MPH, 2-pole motor, No Power Option base with ASAP II Rehab Seat
- Quickie Pulse 6 BC (K0848): 6 MPH, 4-pole motor, No Power Option base with ASAP II Rehab Seat
- Quickie Pulse 5 CC (K0849): 5 MPH, 2-pole motor, No Power Option base with Captain's Seat
- Quickie Pulse 6 CC (K0849): 6 MPH, 4-pole motor, No Power Option base with Captain's Seat
- Quickie Pulse 6 SC (K0856): 6 MPH, 4-pole motor, Single Power Option Tilt (SPOT) or Power Recline with ASAP II Rehab Seat or 3.7 Recline Seat
- Quickie Pulse 6 MPC (K0861): 6 MPH, 4-pole motor, Single Power Option Tilt (SPOT) and Power Recline with 3.7 Recline Seat



## Compact Power Base

- Intuitive mid-wheel drive handling
- Small footprint and compact design for indoor accessibility
- Flexible, independent suspension system for moderate outdoor performance
- Linix 2- or 4-pole motors



## Seating Options

### Captain's Seat (CC)

- Convertible seat pan for flexibility in seating options
- Semi-recline and fold-down backrest
- Seat size, legrest, and armrest adjustability
- Mounting capabilities for O2, crutch, or walker holders

### 3.7 Recline Rehab Seat (SC & MPC)

- Option of Recline Curved back or Integrated JAY J3 Back for seamless functionality with power recline and shear reduction system
- Seat-to-floor heights as low as 16.3"
- Cantilever Reclining Armrests stay parallel to ground when reclined



### ASAP II (BC & SC)

- In-field width, depth, armrest, legrest, and back height adjustability
- Manual reclining backrest for increased function and sitting tolerance
- Optional fold down backrest for easy transport

## Power Options

- 3.7 Power Recline with Power Shear Reduction
- 50° of Single Power Option Tilt (SPOT) to optimize positioning
- 9" Seat Elevate for increased access to cabinets and interaction with peers
- Power Elevating Legrests with Articulation



Phone Mount

Tablet Mount

ReCharge™

## NEW! Options for Consumer Convenience

- **RAM® Phone and Tablet Mounts** allow mobile phones and/or tablets to be connected and disconnected for easy viewing
- **Active Controls Power Buddy** provides an auxiliary USB power source for charging devices and optional LED lights for improved visibility
- **RAM® Mounts Self-Leveling Cup Holder** uses gimbal system to keep drinks level during travel or when stopped on inclines
- **Active Controls ReCharge™ with Socket** eliminates the challenges of connecting and disconnecting the battery charger by providing an easy-to-reach adapter

## Quickie Pulse Specifications

SPECIFICATIONS	PULSE 5 BC	PULSE 6 BC	PULSE 5 CC	PULSE 6 CC	PULSE 6 SC	PULSE 6 MPC
Drive Wheel Position	Mid	Mid	Mid	Mid	Mid	Mid
HCPCS Code	K0848	K0848	K0849	K0849	K0856	K0861
Speed	5.0 MPH	6.0 MPH	5.0 MPH	6.0 MPH	6.0 MPH	6.0 MPH
Overall Weight	235 lbs.	235 lbs.	250 lbs.	250 lbs.	260 lbs.	345 lbs.
Overall Width	24"	24"	24"	24"	24"	24"
Overall Length	35"	35"	35"	35"	35"	35"
User Weight Capacity	300 lbs.	300 lbs.	300 lbs.	300 lbs.	300 lbs.	300 lbs.
Seat Widths	12" to 20"	12" to 22"	16", 18", 20"	16", 18", 20"	12" to 22"	14" to 22"
Seat Depths	12" to 20"	12" to 22"	15" to 22"	15" to 22"	12" to 22"	16" to 20"
Seat To Floor Height	16" to 20"	16" to 20"	16" to 18" (from pan)	16" to 18" (from pan)	16.3" to 21.6"	16.3" to 18.3"
Seating Frame Options	ASAP II Rehab Seat	ASAP II Rehab Seat	Captain's Seat	Captain's Seat	ASAP II Rehab Seat or 3.7 Recline	3.7 Recline
Transit Approved	ISO	ISO	WC-19*	WC-19*	Tilt: ISO/Recline: WC-19*	WC-19*
Warranty	Limited/Lifetime	Limited/Lifetime	Limited/Lifetime	Limited/Lifetime	Limited/Lifetime	Limited/Lifetime

\* Conforms with ANSI/RESNA WC-19 standards. Please refer to product Owner's Manual for additional information. Owner's Manuals are available for download at [www.sunrisemedical.com](http://www.sunrisemedical.com)

## *Live Without Limits*

The Quickie S-636 features innovative, Group 3 performance and versatile seating and positioning options to accommodate both general and specialized rehabilitation needs.

# QUICKIE® S-6 Series



### **Available in Six Group 3 and 4 Models:**

- Quickie S-636 BC (K0848): No Power Option S-636 base with Rehab seat
- Quickie S-636 SC (K0856): Single Power Option (tilt or recline) S-636 base with Rehab or Recline seat
- Quickie S-636 MPC (K0861): Multi Power Option (tilt and recline) S-636 base with Rehab Recline seat
- Quickie S-646 SE BD (K0868): No Power Option S-646 SE base with Rehab seat
- Quickie S-646 SE SD (K0877): Single Power Option (tilt or recline) S-646 SE base with Rehab or Recline seat
- Quickie S-646 SE MPD (K0884): Multi Power Option (tilt and recline) S-646 SE base with Rehab Recline seat





## Hybrid Base with Dynamic Suspension

- Linix 6.5 mph or 8.5 mph, 4-pole standard motor package for all models
- Short wheel base with hybrid wheel placement for mid-wheel maneuverability and rear-wheel all-terrain performance
- Adjustable dampening, independent wheel suspension that provides a smooth ride, even over obstacles
- PosiDrive™ technology increases the overall strength of the drive train
- Optional LED headlight, blinker and tail light package for safety
- Painted frame available in 23 different color options

## Seating Options



### Rehab Seat

- TRAX System with easy dovetail track adjustments of the seat and back
- Extensive seating and positioning options
- Articulating vent tray option

### Power Recline Seat

- Wide range of seating and positioning options
- Low seat to floor options with and without Power CG Tilt
- Headrest design is easily accessed and does not require tools for removal



## Power Options

- 53° of Power CG Tilt for optimized positioning
- Power Recline with Power Shear Reduction
- Power Elevating Legrests with Articulation

## Quickie S-6 Series Specifications

SPECIFICATIONS	S-636 BC	S-636 SC	S-636 MPC	S-646 SE BD	S-646 SE SD	S-646 SE MPD
HCPCS Code	K0848	K0856	K0861	K0868	K0877	K0884
Drive Wheel Position	RWD	RWD	RWD	RWD	RWD	RWD
Speed	6.5 MPH*	6.5 MPH*	6.5 MPH**	8.5 MPH	8.5 MPH	8.5 MPH
Overall Weight	128 lbs. w/o batteries	128 lbs. w/o batteries	128 lbs. w/o batteries	128 lbs. w/o batteries	128 lbs. w/o batteries	128 lbs. w/o batteries
Overall Width	27.75"	27.75"	27.75"	25.75"	25.75"	25.75"
Overall Length	33.5"	33.5"	33.5"	33.5"	33.5"	33.5"
User Weight Capacity	300 lbs.	300 lbs.	300 lbs.	300 lbs.	300 lbs.	300 lbs.
Seat Widths	14" to 24"	14" to 24"	14" to 24"	14" to 24"	14" to 24"	14" to 24"
Seat Depths	14" to 22"	14" to 22"	14" to 22"	14" to 22"	14" to 22"	14" to 22"
Front Seat To Floor Height	18.5" to 20.5"	18.5" to 20.5"	18.5" to 20.5"	18" to 21"	18.5" to 20.5"	18.5" to 20.5"
Seating Frame Options	Rehab	Rehab or Rehab Recline	Rehab Recline	Rehab	Rehab or Rehab Recline	Rehab Recline
Power Seating Options	N/A	Til or Recline	Tilt and Recline	N/A	Tilt or Recline	Tilt and Recline

\* Optional 5 mph and 8.5 mph motor packages also available

\*\* Optional 8.5 mph motor package also available

## *Live Without Limits*

The powerful and compact Quickie P-222 SE is fast, simple, and extremely maneuverable indoors and out and features a quick-release power pack and optional manual conversion kit for convenience and transportability.

# QUICKIE® P-222 SE



### **Available in Two Group 4 Models:**

- Quickie P-222 SE BD (K0868): No power option P-222 SE base with Rehab seat
- Quickie P-222 SE SD (K0877): Single power option (tilt) P-222 SE base with Rehab seat



## Fast and Maneuverable Base

- Linix 8.5 mph, 4-pole standard motor package
- PosiDrive™ technology increases the overall strength of the drive train and features easy-to-service hubs
- 26" turning radius for excellent maneuverability
- 4" curb climbing capability



## Seating and Power Options

### Rehab Seat

- Simple seating and positioning options mix
- Swing-away and fixed hanger options
- Optional 53° of CG tilt to optimize positioning



## Style Options

- Optional manual wheelchair conversion kit for user self propulsion
- Optional LED light package includes headlights, blinkers and tail light for safety
- Painted frame available in 23 different color options

## Quickie P-222 SE Specifications

SPECIFICATIONS	P-222 SE BD	P-222 SE SD
HCPCS Code	K0868	K0877
Drive Wheel Position	RWD	RWD
Speed	8.5 MPH	8.5 MPH
Overall Weight	115 lbs. w/o batteries	115 lbs. w/o batteries
Overall Width	24.5"	24.5"
Overall Length	31.5" w/o footplates	31.5" w/o footplates
User Weight Capacity	350 lbs.	350 lbs.
Seat Widths	14" to 24"	14" to 24"
Seat Depths	14" to 20"	14" to 20"
Front Seat To Floor Height	18" to 20"	17.5" to 20.5"
Seating Frame Options	Rehab	Rehab
Power Seating Options	N/A	Tilt

# PG Drives Technology

*Option availability varies by chair and model. Please see order form for detailed configuration information.*



Switch-It Head Array

## VR2 Non-Expandable Controls

- Industry standard non-expandable electronics
- Dependable, highly flexible, and simple to use
- 4-way toggle switch for operation of tilt and power legrests without the need for expandable controls
- Optional Attendant Control for operation by caregiver

## R-net Expandable and Non-Expandable Controls


- 8 uniquely programmable profiles for easy transitions between environments
- Optional assignable buttons for convenient control of common functions
- Versatile joystick options
- Optional LED joystick
- Optional joystick w/color display
  - Color Screen for clear operation of functions
  - On-board programming capabilities
  - Available with mono and stereo jacks for adding external actuator and ON/OFF switch control
  - Available with Toggle, Speed Potentiometer, and stereo jacks for quad-friendly controls
- Expandable OMNI Input Control Module
  - Infra Red (IR) functions for operation of IR controlled devices such as TVs and DVD players
  - Bluetooth Mouse-Mover for PC mouse control
  - On-board programming capabilities
- Full specialty control and switch input capabilities
- Optional Attendant Control for operation by caregiver

## Back Upholstery Options



### QUICKIE® Tension Adjustable Upholstery

- Available with Ballistic Nylon or 3DX Microclimatic fabric
- 1" width adjustment

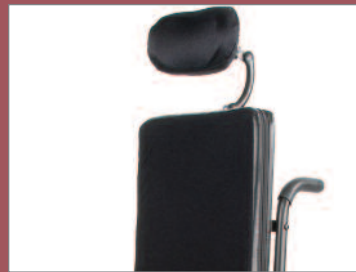

### SPOT Curved Back

- Solid Support
- Optional preassembled Whitmyer headrests
- Optional AES positioning components



### JAY® J3™ Backrest

- Two-point or four-point quick-release hardware
- Range of back heights and contours to meet user needs
- Can mount Whitmyer headrests or lateral positioning supports



### SPO Backrest

- JAY or ART SPO Backrests are designed for the unique needs of a client in power tilt.
- Provides postural management in a package that can be customized in the field.
- Shape, height, width and thickness of the foam customizable by user.
- ART SPO Back offers custom shapes: i-back, scapular cut-out and custom rectangular.

# Electronics

*Option availability varies by chair and model. Please see order form for detailed configuration information.*

## Drive Controls



VR2 Non-Expandable Controller



R-net Expandable Controller



R-net 120amp Non-Expandable/Expandable Controller



R-net Environmental Control Output Module (ECM)



R-net Bluetooth® Mouseover

## Joysticks



Non-expandable VR-2 Joystick



R-net Expandable LED Joystick



Joystick w/Stereo Jacks & Color Display



Joystick w/Toggle, Speed Potentiometer, Stereo Jacks & Color Display



R-net OMNI Input Control Module

## Joystick Mounts



Fixed



Swing-Away




VR-2




R-net

## Attendant Controls

## Joystick Handles

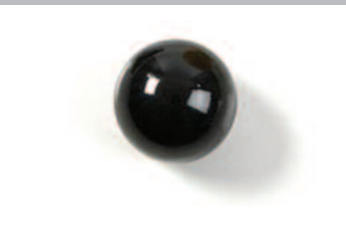

Ball

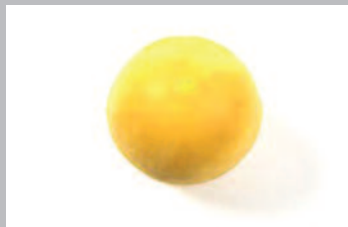

Foam Ball



T-Handle



Goalpost

Option availability varies by chair and model. Please see order form for detailed configuration information.

## Switch Inputs for Joysticks or Omni



Ribbon Switch



Buddy Button Switch



Disc Switch



Wobble Switch



Micro Light Switch



Egg Switch - Black



Stereo to 2 Mono Splitter




Switch-It® Mini 1 Button on Mono Plug



Switch-It® Mini 2 Button on 2 Mono Plugs



Switch-It® Mini 2 Button on 1 Stereo Plug

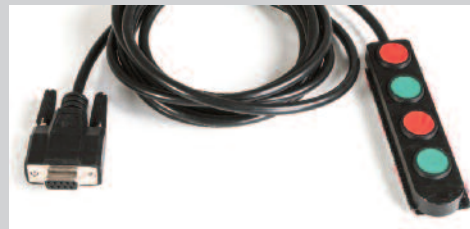

Switch-It® Mini 4 Button on DB9 Plug

## Proportional Controls - Chin Control



Switch-It® Micro-Pilot Bib Chin Control Pkg, Standard 6" or Large 8"



Switch-It® Micro-Guide Bib Chin Control Pkg, Standard 6" or Large 8"



Quickie® Proportional Chin Control Package



ASL® Proportional Micro Extremity Joystick Chin/Lip Control Package

# Proportional Controls - Hand Control



ASL® Proportional Micro Extremity Joystick Midline with Gatlin Package



ASL® Mushroom Joystick Midline with Gatlin Package



ASL® Compact Joystick with Integrated Buttons Midline with Gatlin Package



ASL® Micro Mini Joystick with Hand Pad Midline with Gatlin Package



Switch-It® Micro-Pilot Joystick with Bullet Tray Package



Switch-It® Micro-Guide Joystick with Bullet Tray Package

# Switch Control



ASL® Proximity Head Array Package 3-Switch (Adult or Pediatric)



Switch-It® Head Array Package 3-Switch Curved (Long or Stubby)




Switch-It® Head Array 5-Switch Flat



Quickie® Sip 'N' Puff System

# Multiple Switches

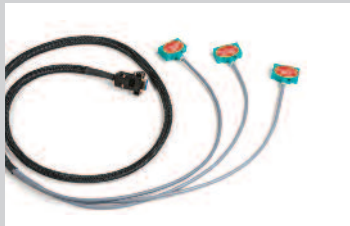

Switch-It® 3 Naked Proximity Switches on DB9 Connector



Switch-It® 4 Naked Proximity Switches on DB9 Connector



Switch-It® 5 Naked Proximity Switches on DB9 Connector



ASL® 5-Switch Adaptor Box on DB9 Connector

# Programmers



VR-2 Hand Held Programmer



R-net PC Programmer



DTT VR2 and R-net

# Base Options


*Option availability varies by chair and model. Please see order form for detailed configuration information.*



**Knobby Drive Wheels**



**4-Point Transit Tie-Downs**

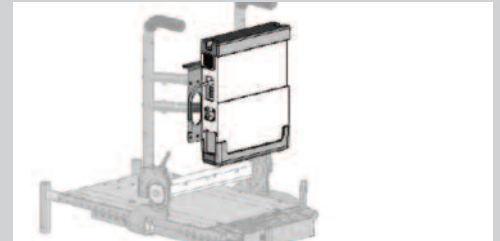

**Fixed Vent Tray**



**Gp 24 Batteries with Std Posts**



**Gp 24 Batteries with Flat Tops**



**22NF Batteries with Std Posts**



**22NF Batteries with Flat Tops**



**Gp 34 Batteries with Flat Tops**



**8A Dual Mode (Off-Board) CTE**




# Accessories

*Option availability varies by chair and model. Please see order form for detailed configuration information.*



RAM® Cup Holder

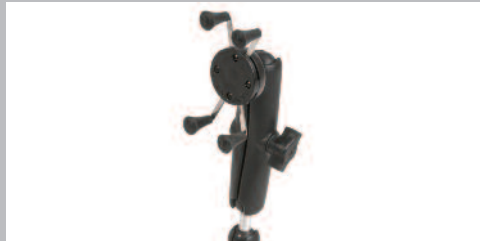

RAM® Phone Mount



RAM® Tablet Mount



Active Controls Power Buddy



Active Controls ReCharge



Switch-It® Remote Stop Switch



Auto Style Buckle



Belt 2" Aircraft Buckle



Belt 2" Aircraft Buckle Padded



O2 Holder



Crutch Holder



Walker Holder



Tuffak Tray Table



Black Utility Backpack



Black/Charcoal Kid's Backpack

# Armrests

*Option availability varies by chair and model. Please see order form for detailed configuration information.*



Dual Post Flip Back  
Height Adjustable



Flip Back



Standard Cantilever with 2 Point  
Locking Post



Standard Cantilever



Low Cantilever



Standard Single Post



Reclining Standard Height



Cantilever Reclining

## Armpads



Waterfall Armpad



Standard Armpad



Reclining Armpad



Ottobock® Reclining Armpad



Ottobock® Flat Hand Pad

## Positioning Options



\*Whitmyer® Plush  
with Lynx hardware



\*Whitmyer® Plush  
with Cobra hardware



\*Whitmyer® Plush  
with Axis hardware



Chest Strap




Lateral Thoracic Supports



JAY® J3™ Adjustable Lateral



Swing-Away Knee Adductor

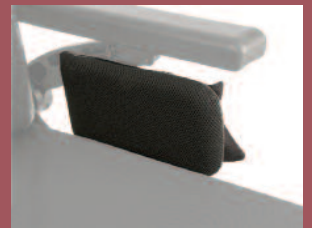

Long Thigh Support  
with mesh pocket

*\*Various Whitmyer head supports are available with the hardware above.*

# Legrests

*Option availability varies by chair and model. Please see order form for detailed configuration information.*



65° & 70° Swing-Away



75° Tapered Swing-Away



90° Swing-Away



Manual ELR Swing-Away



Fixed Center Mount



Centermount Power ELR  
with Articulation



Swing Away Power ELR  
with Articulation



Kids Angle Adjust with Leg Strap

## Footplates




Composite with Heel Loop



Adult Angle Adjust with  
Heel Loop



Kids Angle Adjust with Leg Strap

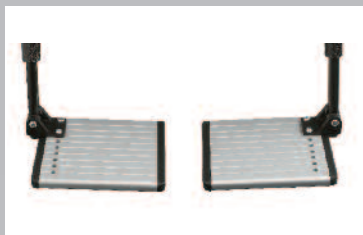

Aluminum



One Piece Angle Adjust Footplate



Pediatric Footplate



Centermount Single/Dual Footplate



www.built4me.us

## Product Customizations


Our Built-4-Me Team is dedicated to creating product customizations beyond what is listed on our product order forms. Experts in the engineering and manufacturing of custom rehab products will listen to your specific requests and build a chair as individual as your fingerprint!

### Featured Power Chair Customizations


- 24" Seat width on QM-715 HD
- Captain's Seat on QM-710

# JAY<sup>®</sup> Seating Selection


## JAY CUSHIONS

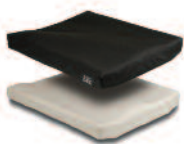


**Basic**  
E2601, E2602

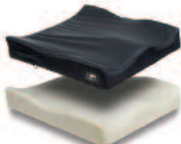

**BasicPRO™**  
E2601, E2602

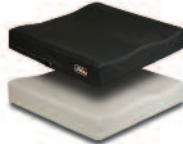

**GO**  
E2603, E2604

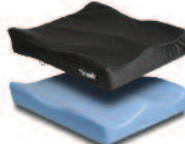


**Soft Combi™ P**  
E2605, E2606



**Care**  
E2607, E2608



**Zip™**  
E2607

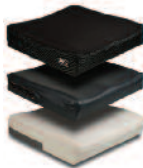

**Ion™**  
E2607, E2608



**Easy**  
E2607, E2608



**Union™**  
E2307, E2608



**X2™**  
E2622



**J2™**  
E2622, E2623,  
E2624, E2625



**Fusion™**  
E2622, E2623



**J3™**  
E2622, E2623,  
E2624, E2625

## JAY BACKS



**Basic**  
E2611



**GO**  
E2611, E2612



**Care**  
E2615



**Focus Point**  
E2620



**Zip™**  
E2613



**J3™ Posterior**  
E2613, E2614



**J3™ Posterior & Lateral**  
E2615



**J3™ Posterior & Deep Lat/Cont**  
E2620